

PAGANA • PAGANA • PAGANA

MOSBY'S  
DIAGNOSTIC  
& LABORATORY  
TEST REFERENCE

---

Fourteenth Edition

ELSEVIER

## ROUTINE BLOOD TESTING

*Many laboratory tests include the direction to perform routine blood testing. The protocol for those tests is presented here and is cross-referenced within the many tests requiring them.*

### Before

- Follow proper patient identification protocols to avoid wrong patient events. Usually name and date of birth are used as two identifiers.

**PT** Explain the procedure to the patient.

**PT** Tell the patient if fasting is necessary. (Fasting is most commonly required with glucose and lipid studies.)

**PT** If fasting is required, instruct the patient not to consume any food or fluids. Only water is permitted. Fasting requirements usually vary from 8 to 12 hours.

**PT** Instruct the patient to continue taking medications unless told otherwise by the health-care provider.

### During

- Collect the blood in a properly color-coded test tube (Table A, p. xiv), which indicates the presence or absence of additives. Tube stopper colors may vary with different manufacturers. If uncertain, verify with the laboratory.

### After

- Apply pressure or a pressure dressing to the venipuncture site.
- Assess the site for bleeding.

**PT** = Patient teaching

## ROUTINE URINE TESTING

*Many laboratory tests include the direction to perform routine urine testing. The protocol for those tests is presented here and is cross-referenced within the many tests requiring them.*

### **Before**

- Follow proper patient identification protocols to avoid wrong patient events. Usually, name and date of birth are used as two identifiers.

**PT** Explain the procedure to the patient.

**PT** Inform the patient if food or fluid restrictions are needed.

### **During**

#### **Random, fresh, or spot specimen**

**PT** Instruct the patient to urinate into an appropriate nonsterile container.

#### **24-hour specimen**

1. Begin the 24-hour collection by discarding the first specimen.
2. Collect all urine voided during the next 24 hours.
3. Show the patient where to store the urine.
4. Keep the urine on ice or refrigerated during the collection period. Foley bags are kept in a basin of ice. Some collections require a preservative. Check with the laboratory.
5. Post the hours for the urine collection in a prominent place to prevent accidentally discarding a specimen.
6. Instruct the patient to void before defecating so that urine is not contaminated by stool.
7. Remind the patient not to put toilet paper in the urine collection container.
8. Collect the last specimen as close as possible to the end of the 24-hour period. Add this urine to the collection.

### **After**

- Transport the specimen promptly to the laboratory.

**PT** = Patient teaching

## COMMON REFERENCE RANGES FOR HEALTHY ADULTS

*(NOTE: These values are generalizations. Each laboratory has specific ranges.)*

<b>Blood Count/Hematology</b>	<b>Page</b>
WBC: $5-10 \times 10^9/L$	974
RBC: ♂ $4.7 - 6.1 \times 10^6/\mu L$ ; ♀ $4.2-5.4 \times 10^6/\mu L$	770
Hemoglobin: ♂ 14-18 g/dL; ♀ 12-16 g/dL	488
Hematocrit: ♂ 40%-52%; ♀ 36%-47%	485
Platelets: $150-400 \times 10^9/L$	706
Prothrombin time (PT): 11-12.5 sec	753
International normalized ratio (INR): 0.8-1.1	753
Activated partial thromboplastin time (APTT): 30-40 sec	681
Glycosylated hemoglobin: 4%-5.9%	471
 <b>Electrolytes and Gastrointestinal, Renal, and Liver Function</b>	
Sodium: 136-145 mEq/L	835
Potassium: 3.5-5 mEq/L	724
Chloride: 98-106 mEq/L	233
CO <sub>2</sub> content (bicarbonate): 23-30 mEq/L	197
Blood urea nitrogen (BUN): 10-20 mg/dL	155
Creatinine: 0.5-1.1 mg/dL	301
Glucose: 74-106 mg/dL	462
Calcium: 9-10.5 mg/dL	189
Amylase: 60-120 U/L	56
Lipase: 0-160 U/L	562
Protein (total): 6.4-8.3 g/dL	476
Albumin: 3.5-5 g/dL	476
Bilirubin (total): 0.3-1 mg/dL	137
Bilirubin (direct): 0.1-0.3 mg/dL	137
Alkaline phosphatase (ALP): 30-120 U/L	29
Alanine aminotransferase (ALT): 4-36 U/L	21
Aspartate aminotransferase (AST): 0-35 U/L	125
Gamma-glutamyl transpeptidase (GGT): 8-38 U/L	435
 <b>Lipids</b>	
Triglycerides: 40-180 mg/dL	908
Total cholesterol: < 200 mg/dL	235
High-density lipoproteins (HDL): ♂ > 45 mg/dL; ♀ > 55 mg/dL	565
Low-density lipoproteins (LDL): < 130 mg/dL	565



**MOSBY'S**  
**DIAGNOSTIC AND  
LABORATORY  
TEST REFERENCE**

---

*Fourteenth Edition*

**Kathleen Deska Pagana, PhD, RN**

Professor Emeritus  
Department of Nursing  
Lycoming College  
Williamsport, Pennsylvania  
President, Pagana Keynotes & Presentations  
<http://www.KathleenPagana.com>

**Timothy J. Pagana, MD, FACS**

Medical Director Emeritus  
The Kathryn Candor Lundy Breast Health Center  
Susquehanna Health System  
Williamsport, Pennsylvania

**Theresa Noel Pagana, MD, FAAEM**

Emergency Medicine Physician  
Virtua Voorhees Hospital  
Voorhees, New Jersey

**ELSEVIER**

# ELSEVIER

3251 Riverport Lane  
St. Louis, Missouri 63043

MOSBY'S DIAGNOSTIC AND LABORATORY  
TEST REFERENCE, FOURTEENTH EDITION ISBN: 978-0-323-60969-2

Copyright © 2019 by Elsevier, Inc. All rights reserved.

Previous editions copyrighted 2017, 2015, 2013, 2011, 2009, 2007, 2005, 2003,  
2001, 1999, 1997, 1995, 1992

No part of this publication may be reproduced or transmitted in any form or by any means, electronic or mechanical, including photocopying, recording, or any information storage and retrieval system, without permission in writing from the publisher. Details on how to seek permission, further information about the Publisher's permissions policies and our arrangements with organizations such as the Copyright Clearance Center and the Copyright Licensing Agency, can be found at our website: [www.elsevier.com/permissions](http://www.elsevier.com/permissions).

This book and the individual contributions contained in it are protected under copyright by the Publisher (other than as may be noted herein).

## Notices

Practitioners and researchers must always rely on their own experience and knowledge in evaluating and using any information, methods, compounds or experiments described herein. Because of rapid advances in the medical sciences, in particular, independent verification of diagnoses and drug dosages should be made. To the fullest extent of the law, no responsibility is assumed by Elsevier, authors, editors or contributors for any injury and/or damage to persons or property as a matter of products liability, negligence or otherwise, or from any use or operation of any methods, products, instructions, or ideas contained in the material herein.

Library of Congress Control Number: 2018952399

*Senior Content Strategist:* Yvonne Alexopoulos

*Content Development Manager:* Lisa P. Newton

*Senior Content Development Specialist:* Tina Kaemmerer

*Publishing Services Manager:* Deepthi Unni

*Project Manager:* Janish Ashwin Paul

*Design Direction:* Renee Duenow

Printed in United States

Last digit is the print number: 9 8 7 6 5 4 3 2 1



Working together  
to grow libraries in  
developing countries

[www.elsevier.com](http://www.elsevier.com) • [www.bookaid.org](http://www.bookaid.org)

*With love and adoration,  
we dedicate this book to our grandchildren:*

Ella Marie Gaul  
Jocelyn Elizabeth Gaul  
Timothy William Gaul  
Justin Aquinas Gaul  
Juliana Kathleen Pericci  
Luke Michael Pericci  
John Henry Bullen V  
Hunter Timothy Bullen

—KDP, TJP

*Thanks for your love and support,  
John Bullen IV*

—TNP

## reviewers

### **Brenda Barnes, BS, MS, PhD**

Director, Medical Lab Science  
Program;  
Director, EdD Health  
Professions Education  
Program;  
Professor  
Allen College-UnityPoint  
Health  
Waterloo, Iowa

### **Tammy R. Dean, RN, BSN**

Program Director  
The Prince William County  
School of Practical Nursing  
Manassas, Virginia

### **Sue Ellen Edrington, MSN, RN**

Assistant Professor  
Leighton School of Nursing  
Marian University  
Indianapolis, Indiana

### **Lorraine Kelley, DNP, RN**

Nursing Instructor  
Pensacola State College  
Cantonment, Florida

### **Marilyn Kelly, RN, LNC**

Health Sciences Instructor  
Discovery Community  
College  
Nanaimo, British Columbia,  
Canada

### **Peter Miskin, RN, DHSc**

Adjunct Assistant Professor  
Samuel Merritt University,  
Peninsula Learning Center  
San Mateo, California

### **Flora Sayson**

Instructor, Nursing  
College of Southern Nevada  
Las Vegas, Nevada

### **Jessica Massengill, MHA-Edu, BSN, RN**

Coordinator of Health  
Sciences  
Tennessee College of Applied  
Technology-Harriman  
Harriman, Tennessee

### **Shopha Tserotas, MS, RN**

Coordinator, Weekend  
Program  
Assistant Clinical Professor  
Texas Woman's University  
Dallas, Texas

## preface

The 14th edition of *Mosby's Diagnostic and Laboratory Test Reference* provides the user with an up-to-date, essential reference that allows easy access to clinically relevant laboratory and diagnostic tests. A unique feature of this handbook is its consistent format, which allows for quick reference without sacrificing the depth of detail necessary for a thorough understanding of diagnostic and laboratory testing. All tests begin on a new page and are listed in alphabetical order by their complete names. The alphabetical format is a strong feature of the book; it allows the user to locate tests quickly without first having to place them in an appropriate category or body system. The User's Guide to Test Preparation and Procedures section outlines the responsibilities of health-care providers to ensure that the tests are accurately and safely performed. Use of this guide should eliminate the need for test repetition resulting from problems with patient preparation, test procedures, or collection techniques. Information on radiation exposure and risks has been added. Every feature of this book is designed to provide pertinent information in a sequence that best simulates priorities in the clinical setting.

The following information is provided, wherever applicable, for effective diagnostic and laboratory testing:

- **Name of test.** Tests are listed by their complete names. A complete list of abbreviations and alternate test names follows each main entry.
- **Type of test.** This section identifies whether the test is, for example, an x-ray procedure, ultrasound, nuclear scan, blood test, urine test, sputum test, or microscopic examination of tissue. This section helps the reader identify the source of the laboratory specimen or location of the diagnostic procedure.
- **Normal findings.** Where applicable, normal values are listed for the infant, child, adult, and elderly person. Also, where appropriate, values are separated into male and female. It is important to realize that normal ranges of laboratory tests vary from institution to institution. This variability is even more obvious among the various laboratory textbooks. For this reason, we have deliberately chosen not to add a table of normal values as an appendix, and we encourage the user to check the normal values at the institution where the test is performed. This should be relatively easy because laboratory reports include normal values. Results are given in both conventional units and the International System of Units (SI units) where possible.

## vi preface

- **Possible critical values.** These values give an indication of results that are well outside the normal range. These results require health-care provider notification and usually result in some type of intervention. The Joint Commission is looking at the timely and reliable communication of critical laboratory values as one of its patient safety goals.
- **Test explanation and related physiology.** This section provides a concise yet comprehensive description of each test. It includes fundamental information about the test itself, specific indications for the test, how the test is performed, what disease or disorder the various results may show, how it will affect the patient or client, and relevant pathophysiology that will enhance understanding of the test.
- **Contraindications.** These data are crucial because they alert health-care providers to patients to whom the test should not be administered. Patients highlighted in this section frequently include those who are pregnant, are allergic to iodinated or contrast dyes, or have bleeding disorders.
- **Potential complications.** This section alerts the user to potential problems that necessitate astute assessments and interventions. For example, if a potential complication is renal failure, the implication may be to hydrate the patient before the test and force fluids after the test. A typical potential complication for many x-ray procedures is allergy to iodinated dye. Patient symptoms and appropriate interventions are described in detail.
- **Interfering factors.** This section contains pertinent information because many factors can invalidate the test or make the test results unreliable. An important feature is the inclusion of drugs that can interfere with test results. Drugs that increase or decrease test values are always listed at the end of this section for consistency and quick access. A drug symbol (☒) is used to emphasize these drug interferences.
- **Procedure and patient care.** This section emphasizes the role of nurses and other health-care providers in diagnostic and laboratory testing by addressing psychosocial and physiologic interventions. Patient teaching priorities are noted with a special icon (PT) to highlight information to be communicated to patients. For quick access to essential information, this section is divided into before, during, and after time sequences.
  - **Before.** This section addresses the need to explain the procedure and to allay patient concerns or anxieties. If patient consent is usually required, this is listed as a bulleted item.

Other important features include requirements such as fasting, obtaining baseline values, and performing bowel preparations. Radiation risk is addressed with x-rays and nuclear medicine studies.

- **During.** This section gives specific directions for clinical specimen studies (e.g., urine and blood studies). Diagnostic procedures and their variations are described in a numbered, usually in a step-by-step format. Important information, such as who performs the test, where the study is performed, patient sensation, and duration of the procedure, is bulleted for emphasis. The duration of the procedure is very helpful for patient teaching because it indicates the time generally allotted for each study.
- **After.** This section includes vital information that the nurse or other health-care provider should heed or convey after the test. Examples include such factors as maintaining bed rest, comparing pulses with baseline values, encouraging fluid intake, and observing the patient for signs and symptoms of sepsis.
- **Abnormal findings.** As the name implies, this section lists the abnormal findings for each study. Diseases or conditions that may be indicated by increased (▲) or decreased (▼) values are listed where appropriate.
- **Notes.** This blank space at the end of the tests facilitates individualizing the studies according to the institution at which the test is performed. Variations in any area of the test (e.g., patient preparation, test procedure, normal values, postprocedural care) can be noted here.

This logical format emphasizes clinically relevant information. The clarity of this format allows for quick understanding of content essential to both students and health-care providers. Color has been used to help locate tests and to highlight critical information (e.g., possible critical values). Color is also used in the illustrations to enhance the reader's understanding of many diagnostic procedures (e.g., bronchoscopy, fetoscopy, endoscopic retrograde cholangiopancreatography [ERCP], pericardiocentesis, transesophageal echocardiography [TEE]). Many tables are used to simplify complex material on such topics as bioterrorism infectious agents, blood collection tubes, hepatitis testing, and protein electrophoresis. Extensive cross-referencing exists throughout the book, which facilitates understanding and helps the user tie together or locate related studies, such as hemoglobin and hematocrit.



## viii preface

Standard guidelines for routine blood and urine testing are located on the inside front cover for easy access. A list of abbreviations for test names is included on the book's endpapers.

**Appendix A** includes a list of studies according to *body system*. This list may familiarize the user with other related studies the patient or client may need or the user may want to review. This should be especially useful for students and health-care providers working in specialized areas.

**Appendix B** provides a list of studies according to *test type*. This list may help the user read and learn about similarly performed tests and procedures (e.g., barium enema and barium swallow).

**Appendix C** provides a list of blood tests used for *disease and organ panels*.

**Appendix D** provides a list of *symbols and units of measurement*.

Finally, a comprehensive index includes the names of all tests, their synonyms and abbreviations, and any other relevant terms found in the tests.

New to this edition is a table of Common Reference Ranges added to the inside front cover. This adds to the user-friendly aspect of this book by quickly identifying common reference ranges. This is a good starting point for students and a quick reference for routine lab values. However, because lab values vary from institution to institution, be sure to use the normal values of the lab performing the test.

Many new studies, such as alpha defensin, ceramides, and small intestinal bacterial overgrowth tests, have been added. All other studies have been revised and updated. Outdated studies have been eliminated.

We sincerely thank our editors for their enthusiasm and continued support. We are most grateful to the many nurses and other health-care providers who made the first 13 editions of this book so successful. Thank you so much. This success validated the need for a user-friendly and quick-reference approach to laboratory and diagnostic testing.

We sincerely invite additional comments from current users of this book so that we may continue to provide useful, relevant diagnostic and laboratory test information to users of future editions.

**Kathleen D. Pagana**  
**Timothy J. Pagana**  
**Theresa N. Pagana**

## contents

**Routine blood testing, inside front cover**

**Routine urine testing, inside front cover**

**Common reference ranges, inside front cover**

**List of figures, x**

**User's guide to test preparation and procedures, xi**

**Diagnostic and laboratory tests, 1**

*Tests presented in alphabetical order*

### **Appendices**

Appendix A: List of tests by body system, 988

Appendix B: List of tests by type, 1000

Appendix C: Disease and organ panels, 1011

Appendix D: Symbols and units of measurement, 1015

**Bibliography, 1017**

**Index, 1019**

**Abbreviations for diagnostic and laboratory tests, inside back cover**

## list of figures

- Figure 1** Ultrasound of the abdomen, 2
- Figure 2** Amniocentesis, 53
- Figure 3** Immunofluorescent staining of antinuclear antibodies, 88
- Figure 4** Arthroscopy, 122
- Figure 5** Bilirubin metabolism and excretion, 138
- Figure 6** Bone marrow aspiration, 163
- Figure 7** Bronchoscopy, 185
- Figure 8** Cardiac catheterization, 204
- Figure 9** Chorionic villus sampling, 242
- Figure 10** Hemostasis and fibrinolysis, 251
- Figure 11** Colposcopy, 262
- Figure 12** Cardiac enzymes after myocardial infarction, 298
- Figure 13** Cystoscopic examination of the male bladder, 309
- Figure 14** Ureteral catheterization through the cystoscope, 310
- Figure 15** Disseminated intravascular coagulation, 331
- Figure 16** Ductoscopy, 337
- Figure 17** ECG planes of reference, 343
- Figure 18** Electrocardiography, 344
- Figure 19** Endoscopic retrograde cholangiopancreatography, 366
- Figure 20** Esophageal function studies, 380
- Figure 21** Fetoscopy, 420
- Figure 22** Glucose tolerance test, 468
- Figure 23** Hematocrit, 486
- Figure 24** Holter monitoring, 510
- Figure 25** Hysteroscopy, 533
- Figure 26** Liver biopsy, 571
- Figure 27** Lumbar puncture, 581
- Figure 28** Transbronchial needle biopsy, 585
- Figure 29** Stereotactic breast biopsy, 607
- Figure 30** Oximetry, 659
- Figure 31** Papanicolaou (Pap) smear, 669
- Figure 32** Paracentesis, 672
- Figure 33** Pericardiocentesis, 690
- Figure 34** Rectal ultrasonography, 741
- Figure 35** Lung volumes and capacities, 760
- Figure 36** Renal biopsy, 777
- Figure 37** Renovascular hypertension, 784
- Figure 38** Rectal culture of the female, 816
- Figure 39** Urethral culture of the male, 817
- Figure 40** Thoracentesis, 865
- Figure 41** Fibrin clot formation, 875
- Figure 42** Transesophageal echocardiography, 904

## user's guide to test preparation and procedures

Health-care economics demands that laboratory and diagnostic testing be performed accurately and in the least amount of time possible. Tests should not have to be repeated because of improper patient preparation, test procedure, or specimen collection technique. Patient identification protocols should be followed to avoid wrong patient events. Two patient identifiers, such as name and date of birth, are usually used. The following guidelines delineate the responsibilities of health-care providers to ensure safety of test procedures and accuracy of test results. Guidelines are described for the following major types of tests: blood, urine, stool, x-ray, nuclear scanning, ultrasound, and endoscopy.

### Blood tests

#### Overview

Blood studies are used to assess a multitude of body processes and disorders. Common studies include enzymes, serum lipids, electrolyte levels, red and white blood cell counts, clotting factors, hormone levels, and levels of breakdown products (e.g., blood urea nitrogen).

*Multiphasic screening machines* can perform many blood tests simultaneously using a very small blood sample. The advantages of using these machines are that results are available quickly and the cost is lower when compared with individually performing each test.

Appendix C provides a list of current disease and organ panels. For example, the basic metabolic panel and the comprehensive metabolic panel have replaced the Chem-7 and Chem-12 panels. These changes are the result of federal guidelines that have standardized the nomenclature for chemistry panels.

#### Guidelines

- Observe universal precautions when collecting a blood specimen.
- Check whether fasting is required. Many studies, such as fasting blood sugar and cholesterol levels, require fasting for a designated period of time. Water is permitted.
- If ordered, withhold medications until the blood is drawn.
- Record the time of day when the blood test is drawn. Some blood test results (e.g., those for cortisol) vary according to a diurnal pattern, and this must be considered when blood levels are interpreted.

## xii user's guide to test preparation and procedures

- In general, two or three blood tests can be done per tube of blood collected (e.g., two or three chemistry tests from one red-top tube of blood).
- Note the patient's position for certain tests (e.g., renin, because levels are affected by body position).
- Collect the blood in a properly color-coded test tube. Blood collection tubes have color-coded stoppers to indicate the presence or absence of different types of additives (preservatives and anticoagulants). A preservative prevents change in the specimen, and an anticoagulant inhibits clot formation or coagulation. Charts are available from the laboratory indicating the type of tube needed for each particular blood test. A representative chart is shown in [Table A](#), p. xiv.
- Follow the recommended *order of draw* when collecting blood in tubes. Draw specimens into nonadditive (e.g., red-top) tubes before drawing them into tubes with additives. This prevents contamination of the blood specimen with additives that may cause incorrect test results. Fill the tubes in the following order:
  1. Blood culture tubes (to maintain sterility)
  2. Nonadditive tubes (e.g., red-top)
  3. Coagulation tubes (e.g., blue-top)
  4. Heparin tubes (e.g., green-top)
  5. Ethylenediaminetetraacetic acid (EDTA) tubes (e.g., lavender-top)
  6. Oxalate/fluoride tubes (e.g., gray-top)
- To obtain valid results, do not fasten the tourniquet for longer than 1 minute. Prolonged tourniquet application can cause stasis and hemoconcentration.
- Collect the blood specimen from the arm without an intravenous (IV) device, if possible. IV infusion can influence test results.
- Do not use the arm bearing a dialysis arteriovenous fistula for venipuncture unless the physician specifically authorizes it.
- Because of the risk of cellulitis, do not take specimens from the side on which a mastectomy or axillary lymph node dissection was performed.
- Follow the unit guidelines for drawing blood from an indwelling venous catheter (e.g., a triple-lumen catheter). Guidelines will specify the amount of blood to be drawn from the catheter and discarded before blood is collected for laboratory studies. The guidelines will also indicate the amount and type of solution needed to flush the catheter after drawing the blood to prevent clotting.

**TABLE A** Common blood collection tubes

Top color	Additive	Purpose	Test examples
Red	Clot activator	Allows blood sample to clot Separates the serum for testing	Chemistry Bilirubin Blood urea nitrogen
Red/black	Clot activator and gel for serum separator	Serum separator tube for serum determinatives in chemistry and serology	Chemistry, serology
Royal blue	Heparin/EDTA	Provides low levels of trace elements	Trace metals, toxicology
Tan	Heparin/EDTA	Contains no lead	Lead determinatives
Purple or lavender	EDTA	Prevents blood from clotting	Hematology CBC
Gray	Oxalate/fluoride	Prevents glycolysis	Chemistry Glucose
Green	Heparin	Prevents blood from clotting when plasma needs to be tested	Chemistry Ammonia
Blue (light)	Sodium citrate	Prevents blood from clotting when plasma needs to be tested	Prothrombin time Partial thromboplastin time
Black	Sodium citrate	Binds calcium to prevent blood clotting	Westergren ESR
Yellow	Citrate dextrose	Preserves red cells	Blood cultures, blood banking studies

*CBC*, Complete blood count; *EDTA*, ethylenediaminetetraacetic acid; *ESR*, erythrocyte sedimentation rate.

## xiv user's guide to test preparation and procedures

- Do not shake the blood specimen. Hemolysis may result from vigorous shaking and can invalidate test results. Use gentle inversions.
- Collect *blood cultures* before the initiation of antibiotic therapy. Blood cultures are often drawn when the patient manifests a fever. Often two or three cultures are taken at 30-minute intervals from different venipuncture sites.
- *Skin punctures* can be used for blood tests on capillary blood. Common puncture sites include the fingertips, earlobes, and heel surfaces. Fingertips are often used for small children, and the heel is the most commonly used site for infants.
- Ensure that the blood tubes are correctly labeled and delivered to the laboratory.
- After the specimen is drawn, apply pressure or a pressure dressing to the venipuncture site. Assess the site for bleeding.
- If the patient fasted before the blood test, reinstitute the appropriate diet.

## Urine tests

### Overview

Urine tests are easy to obtain and provide valuable information about many body system functions (e.g., kidney function, glucose metabolism, and various hormone levels). The ability of the patient to collect specimens appropriately should be assessed to determine the need for assistance.

### Guidelines

- Observe universal precautions in collecting a urine specimen.
- Use the first morning specimen for routine urinalysis because it is more concentrated. To collect a first morning specimen, have the patient void before going to bed and collect the first urine specimen immediately upon rising.
- *Random* urine specimens can be collected at any time. They are usually obtained during daytime hours and without any prior patient preparation.
- If a *culture and sensitivity (C&S)* study is required or if the specimen is likely to be contaminated by vaginal discharge or bleeding, collect a *clean-catch* or *midstream* specimen. This requires meticulous cleansing of the urinary meatus with an antiseptic preparation to reduce contamination of the specimen by external organisms. Then the cleansing agent must be completely removed because it may contaminate the specimen. Obtain the *midstream* collection by doing the following:



1. Have the patient begin to urinate in a bedpan, urinal, or toilet and then stop urinating. (This washes the urine out of the distal urethra.)
  2. Correctly position a sterile urine container and have the patient void 3 to 4 oz of urine into it.
  3. Cap the container.
  4. Allow the patient to finish voiding.
- One-time *composite* urine specimens are collected over a period that may range anywhere from 2 to 24 hours. To collect a timed specimen, instruct the patient to void and discard the first specimen. This is noted as the *start time* of the test. Instruct the patient to save all subsequent urine in a special container for the designated period. Remind the patient to void before defecating so that urine is not contaminated by feces. Also, instruct the patient not to put toilet paper in the collection container. A preservative is usually used in the collection container. At the end of the specified time period, have the patient void and then add this urine to the specimen container, thus completing the collection process.
  - Collection containers for *24-hour urine specimens* should hold 3 to 4L of urine and have tight-fitting lids. They should be labeled with the patient's name, the starting collection date and time, the ending collection date and time, the name of the test, the preservative, and storage requirements during collection.
  - Many urine collections require preservatives to maintain their stability during the collection period. Some specimens are best preserved by being kept on ice or refrigerated.
  - *Urinary catheterization* may be needed for patients who are unable to void. This procedure is not preferred because of patient discomfort and the risk of patient infection.
  - For patients with an *indwelling urinary catheter*, obtain a specimen by aseptically inserting a needleless syringe into the catheter at a drainage port distal to the sleeve leading to the balloon. Aspirate urine and then place it in a sterile urine container. The urine that accumulates in the plastic reservoir bag should never be used for a urine test.
  - Urine specimens from infants and young children are usually collected in a disposable pouch called a *U bag*. This bag has an adhesive backing around the opening to attach to the child's perineum. After the bag is in place, check the child every 15 minutes to see if an adequate specimen has been collected. Remove the specimen as soon as possible after the collection and then label it and transport it to the laboratory.

## Stool tests

### Overview

The examination of feces provides important information that aids in the differential diagnosis of various gastrointestinal disorders. Fecal studies may also be used for microbiologic studies, chemical determinations, and parasitic examinations.

### Guidelines

- Observe universal precautions in collecting a stool specimen.
- Collect stool specimens in a clean container with a fitted lid.
- Do not mix urine and toilet paper with the stool specimen. Both can contaminate the specimen and alter the results.
- Fecal analysis for occult blood, white blood cells, or qualitative fecal fat requires only a small amount of a randomly collected specimen.
- Quantitative tests for daily fecal excretion of a particular substance require a minimum of a 3-day fecal collection. This collection is necessary because the daily excretion of feces does not correlate well with the amount of food ingested by the patient in the same 24-hour period. Refrigerate specimens or keep them on ice during the collection period. Collect stool in a 1-gallon container.
- A small amount of fecal blood that is not visually apparent is termed *occult blood*. Chemical tests using commercially prepared slides are routinely used to detect fecal blood. Numerous commercial slide tests use guaiac as the indicator. These guaiac tests are routinely done on nursing units and in medical offices.
- Consider various factors (e.g., other diagnostic tests and medications) in planning the stool collection. For example, if the patient is scheduled for x-ray studies with barium sulfate, collect the stool specimen first. Various medications (e.g., tetracyclines and antidiarrheal preparations) affect the detection of intestinal parasites.
- Some fecal collections require dietary restrictions before the collection (e.g., tests for occult blood).
- Correctly label and deliver stool specimens to the laboratory within 30 minutes after collection. If you are unable to deliver the specimen within 30 minutes, it may be refrigerated for up to 2 hours.

## X-ray studies

### Overview

Because of the ability of x-rays to penetrate tissues, x-ray studies provide a valuable picture of body structures. X-ray studies can be as simple as a routine chest x-ray image or as complex as

dye-enhanced cardiac catheterization. With the concern about radiation exposure, it is important to realize that the patient may question if the proposed benefits outweigh the risks involved.

### **Radiation dose**

There are several units used to quantify amount of radiation absorbed from diagnostic imaging tests. The gray (Gy) is the measure of the amount of energy absorbed per unit mass. Because different organs in the body absorb radiation differently, the sievert (Sv) is often used instead of the gray. The sievert is the biological effect of 1 gray of radiation on human body tissue. The sievert is more helpful in comparing radiation exposure to different parts of the body. Radiation doses in medical imaging are typically measured in millisieverts (mSv) or 1/1000 of a sievert. On average, each person receives about 3 mSv of radiation yearly from natural background radiation.

The roentgen equivalent in man (rem) is an older unit to quantify the amount of radiation absorbed from x-rays. 1 rem is equivalent to 0.01 sievert.

See chart below for average amounts of radiation for adults associated with diagnostic testing.

### **Risk of radiation**

Radiation exposure can cause damage to DNA. The body usually rapidly repairs this damage. Mistakes in DNA repair can lead to chromosomal or gene abnormalities that may be linked to cancer induction. The likelihood of cancer induction secondary to radiation exposure increases as the amount of radiation exposure increases. A person has a 5% increase in developing cancer over his or her lifetime after radiation exposure of 1 Sv or more. There can be a lag of many years between radiation exposure and cancer diagnosis. The average lag time is about 10 years after exposure.

The cumulative radiation dose from diagnostic imaging is very small and the benefit of proper diagnosis and treatment of disease generally outweighs the risks. However each patient's current situation and history of radiation must be considered to accurately assess cumulative risks and benefits. Diagnostic procedures with higher radiation doses (e.g., computed tomography [CT] scans) should be clearly justified. Appropriateness Criteria published by the American College of Radiology ([acr.org](http://www.acr.org)) is helpful in justification of performance of x-ray imaging.

Special consideration should be given to pregnant women and children before ordering x-ray imaging because the effects of radiation are more profound in fetuses and young children.

## xviii user's guide to test preparation and procedures

If a woman is pregnant, the risks versus benefits must be carefully considered. Certain studies with low radiation in which the focus of radiation is not on the fetus are obviously safer. (Lead-containing shields can reduce x-ray exposure to fetuses.) Imaging using higher dose of radiation should be given only if the risk of not making the diagnosis is greater than the radiation risk.

Radiation risks are most significant in early fetal period and are less significant as the pregnancy progresses.

Patients with high body mass indexes should also be given extra consideration before ordering imaging studies. These patients often require greater radiation doses to penetrate body thickness to create acceptable images. Nuclear medicines studies are not affected in the same way. Although the x-ray exposure needed to produce one fluoroscopic image is low (compared with radiography), high exposures to patients can result from the time that may be encountered in fluoroscopic procedures.

### Radiation Associated with Diagnostic Testing

Common radiology imaging (XR)	Average adult effective dose <sup>a</sup> (mSv)
Abdomen	0.7
Back (lower)	1.5
Back (upper)	1
Barium enema	8
Bone densitometry (DEXA)	0.001
Cervical spine	0.2
Chest	0.1
Dental	0.005-0.01
Extremity (hands, feet, and so on)	0.001
Fluoroscopy	Per minute
Hip	0.7
Hysterosalpingography	2
Intravenous pyelography (IVP)	3
Mammography	0.4
Neck	0.2
Pelvis	0.6
Skull	0.1
Small bowel follow-through	5
Spine (lumbar)	1.5
Spine (thoracic)	1.0
Upper GI series	6

**Common CT imaging**


---

Abdomen and pelvis	10
Brain (head)	2
Chest	7
Chest (low-dose screening)	2
Coronary angiography	15
CT angiography of the chest	15
Neck	3
Sinuses	0.6
Spine	6
Virtual colonoscopy	10

**Nuclear medicine**


---

Bone scan	5
Brain scan	6.9
Cardiac nuclear stress testing	20-40
Gastric emptying scan	0.4
GI bleeding scan	7.8
Liver scan	3.1
Lung scan (ventilation/perfusion)	2
Parathyroid scan	6.7
Renal scan	2.6
Thyroid scan	4.8
Urea breath test	0.003
WBC scan	6.7

**Other**


---

Abdominal angiogram	12
Cardiac catheterization (diagnostic)	7
Coronary angiogram (stent)	15
Endoscopic retrograde cholangiopancreatography (ERCP)	4
Fluoroscopic Barium Swallow	1.5
Head and neck angiogram	5
Positron emission tomography (PET)/CT	25
Pulmonary angiography	5

---

<sup>a</sup>Effective doses are given as an average, and there may be wide variability in dosing depending on particularities of the test in different testing locations.

## xx user's guide to test preparation and procedures

### Guidelines

- Assess the patient for any similar or recent x-ray procedures.
- Evaluate the patient for *allergies to iodine dye*. Carefully consider the following points:
  1. Many types of contrast media are used in radiographic studies. For example, organic iodides and iodized oils are frequently used.
  2. Allergic reactions to iodinated dye may vary from mild flushing, itching, and urticaria to severe, life-threatening anaphylaxis (evidenced by respiratory distress, drop in blood pressure, or shock). In the unusual event of anaphylaxis, the patient is treated with diphenhydramine (Benadryl), steroids, and epinephrine. Oxygen and endotracheal equipment should be on hand for immediate use.
  3. The patient should always be assessed for allergies to iodine dye before it is administered. Inform the radiologist if an allergy to iodinated contrast is suspected. The radiologist may prescribe Benadryl and steroid preparation to be administered before testing. Usually, hypoallergenic nonionic contrast will be administered to allergic patients during the test.
  4. After the x-ray procedure, evaluate the patient for a delayed reaction to dye (e.g., dyspnea, rashes, tachycardia, hives). This usually occurs within 2 to 6 hours after the test. Treat with antihistamines or steroids.
- Assess the patient for any evidence of dehydration or renal disease. Usually BUN and creatinine tests are obtained before administration of iodine-containing IV contrast. Hydration may be required before the administration of iodine.
- Assess the patient for diabetes. People with diabetes are particularly susceptible to renal disease caused by the administration of iodine-containing IV contrast. Patients with diabetes who take metformin (Glucophage) or glyburide (Micronase) are particularly susceptible to lactic acidosis and hypoglycemia. These medications may be discontinued for 1 to 4 days before and 1 to 2 days after the administration of iodine. Check with the x-ray department.
- Women in their childbearing years should have x-ray examinations during menses or within 10 to 14 days after the onset of menses to avoid possible exposure to a fetus.
- Pregnant women should not have x-ray procedures unless the benefits outweigh the risk of damage to the fetus.

- Note whether other x-ray studies are being planned; schedule them in the appropriate sequence. For example, x-ray examinations that do not require contrast should precede examinations that do require contrast. X-ray studies with barium should be scheduled after ultrasonography.
- Note the necessary dietary restrictions. Such studies as barium enema and intravenous pyelogram (IVP) are more accurate if the patient is kept NPO (fasting from food and liquids) for several hours before the test.
- Determine whether bowel preparations are necessary. For example, barium enemas and IVPs require bowel-cleansing regimens.
- Determine whether signed consent forms are required. These are necessary for most invasive x-ray procedures.
- Remove metal objects (e.g., necklaces, watches) because they can hinder visualization of the x-ray field.
- Patient aftercare is determined by the type of x-ray procedure. For example, a patient having a simple chest x-ray study will not require postprocedure care. However, invasive x-ray procedures involving contrast dyes (e.g., cardiac catheterization) require extensive nursing measures to detect potential complications.

## **Nuclear scanning**

### **Overview**

With the administration of a radionuclide and subsequent measurement of the radiation of a particular organ, functional abnormalities of various body areas (e.g., brain, heart, lung, bones) can be detected. Because the half-lives of the radioisotopes are short, only minimal radiation exposure occurs (See p. xx).

### **Guidelines**

- Radiopharmaceuticals concentrate in target organs by various mechanisms. For example, some labeled compounds (e.g., hippuran) are cleared from the blood and excreted by the kidneys. Some phosphate compounds concentrate in the bone and infarcted tissue. Lung function can be studied by imaging the distribution of inhaled gases or aerosols.
- Note whether the patient has had any recent exposure to radionuclides. The previous study could interfere with the interpretation of the current study.
- Note the patient's age and current weight. This information is used to calculate the dose of radioactive substances.



## xxii user's guide to test preparation and procedures

- Nuclear scans are contraindicated in pregnant women and nursing mothers.
- Many scanning procedures do not require special preparation. However, a few have special requirements. For example, for bone scanning, the patient is encouraged to drink several glasses of water between the time of the injection of the isotope and the actual scanning. For some studies, blocking agents may need to be given to prevent other organs from taking up the isotope.
- For most nuclear scans, a small amount of an organ-specific radionuclide is given orally or injected intravenously. After the radioisotope concentrates in the desired area, the area is scanned. The scanning procedure usually takes place in the nuclear medicine department.
- Instruct the patient to lie still during the scanning.
- Usually encourage the patient to drink extra fluids to enhance excretion of the radionuclide after the test is finished.
- Although the amount of radionuclide excreted in the urine is very low, rubber gloves are sometimes recommended if the urine must be handled. Some hospitals may advise the patient to flush the toilet several times after voiding.

## Ultrasound studies

### Overview

In diagnostic ultrasonography, harmless high-frequency sound waves are emitted and penetrate the organ being studied. The sound waves bounce back to the sensor and are electronically converted into a picture of the organ. Ultrasonography is used to assess a variety of body areas, including the pelvis, abdomen, breast, heart, and pregnant uterus.

### Guidelines

- Most ultrasound procedures require little or no preparation. However, the patient having a pelvic sonogram needs a full bladder, and the patient having an ultrasound examination of the gallbladder must be kept NPO before the procedure.
- Ultrasound examinations are usually performed in an ultrasound room; however, they can be performed in the patient unit.
- For ultrasound, a greasy paste is applied to the skin overlying the desired organ. This paste is used to enhance sound transmission and reception because air impedes transmission of sound waves to the body.

- Because of the noninvasive nature of ultrasonography, no special measures are needed after the study except for helping the patient remove the ultrasound paste.
- Ultrasound examinations have no radiation risk.
- Ultrasound examinations can be repeated as many times as necessary without being harmful to the patient. No cumulative effect has been seen.
- Barium has an adverse effect on the quality of abdominal studies. For this reason, schedule ultrasound of the abdomen before barium studies.
- Large amounts of gas in the bowel obstruct visualization of the bowel. This is because bowel gas is a reflector of sound.

## Endoscopy procedures

### Overview

With the help of a lighted, flexible instrument, internal structures of many areas of the body (e.g., stomach, colon, joints, bronchi, urinary system, and biliary tree) can be directly viewed. The specific purpose and procedure should be reviewed with the patient.

### Guidelines

- Preparation for an endoscopic procedure varies according to the internal structure being examined. For example, examination of the stomach (gastroscopy) will require the passage of an instrument through the esophagus and into the stomach. The patient is kept NPO for 8 to 12 hours before the test to prevent gagging, vomiting, and aspiration. For colonoscopy, an instrument is passed through the rectum and into the colon. Therefore, the bowel must be cleansed and free of fecal material to afford proper visualization. Arthroscopic examination of the knee joint is usually done with the patient under general anesthesia, which necessitates routine preoperative care.
- Schedule endoscopic examinations before barium studies.
- Obtain a signed consent for endoscopic procedures.
- Endoscopic procedures are preferably performed by a physician in a specially equipped endoscopy room or in an operating room. However, some kinds can safely be performed at the bedside.
- Air is instilled into the bowel during colon examinations to maintain patency of the bowel lumen and to afford better visualization. This sometimes causes gas pains.
- In addition to visualization of the desired area, special procedures can be performed. Biopsies can be obtained, and bleeding ulcers can be cauterized. Also, knee surgery can be performed during arthroscopy.

## **xxiv** user's guide to test preparation and procedures

- Specific postprocedure interventions are determined by the type of endoscopic examination performed. All procedures have the potential complication of perforation and bleeding. Most procedures use some type of sedation; safety precautions should be observed until the effects of the sedatives have worn off.
- After colonoscopy and similar studies, the patient may complain of rectal discomfort. A warm tub bath may be soothing.
- Usually keep the patient NPO for 2 hours after endoscopic procedures of the upper gastrointestinal system. Be certain that swallow, gag, and cough reflexes are present before permitting fluids or liquids to be ingested orally.

**abdominal ultrasound** (Abdominal sonogram)

**Type of test** Ultrasound

**Normal findings**

Normal abdominal aorta, liver, gallbladder, bile ducts, pancreas, kidneys, ureters, and bladder

**Test explanation and related physiology**

Ultrasonography provides accurate visualization of the abdominal aorta, liver, gallbladder, pancreas, bile ducts, kidneys, ureters, and bladder. Real-time ultrasound provides an accurate picture of the organ being studied (Figure 1).

The *kidney* is ultrasonographically evaluated to diagnose and locate renal cysts, to differentiate renal cysts from solid renal tumors, to demonstrate renal and pelvic calculi, to document hydronephrosis, to guide a percutaneously inserted needle for cyst aspiration or biopsy, and to place a nephrostomy tube. Ultrasound of the urologic tract is also used to detect malformed or ectopic kidneys and perinephric abscesses. Renal transplantation surveillance is possible with ultrasound.

*Endourethral urologic ultrasound* can also be performed through a stent that has a transducer at its end. The stent probe can be advanced into the bladder where the depth of a tumor into the bladder wall can be measured. In the ureter, stones, tumors, or extraurethral compression can be identified and localized. Finally, in the proximal ureter, renal tumors or cysts can be delineated.


One of the most common uses of ultrasound is the measurement of post void urinary bladder residual. This is a measurement of the amount of urine after micturition. This test can be easily performed at the bedside or in doctor's office with a portable ultrasound unit.

Pelvic, obstetric, prostate, and testes ultrasound are discussed on pp. 685 and 799.

The *abdominal aorta* can be assessed for aneurysmal dilation.

Ultrasound is used to detect cystic structures of the *liver* (e.g., benign cysts, hepatic abscesses, and dilated hepatic ducts) and solid intrahepatic tumors (primary and metastatic). Hepatic ultrasound also can be performed intraoperatively to provide the locations of small, nonpalpable hepatic tumors or abscesses. The *gallbladder* and *bile ducts* can be visualized and examined for evidence of gallstones, polyps, or bile duct dilation. The *pancreas* is examined for evidence of tumors, pseudocysts, acute inflammation, chronic inflammation, or pancreatic abscesses. Because this

## 2 abdominal ultrasound

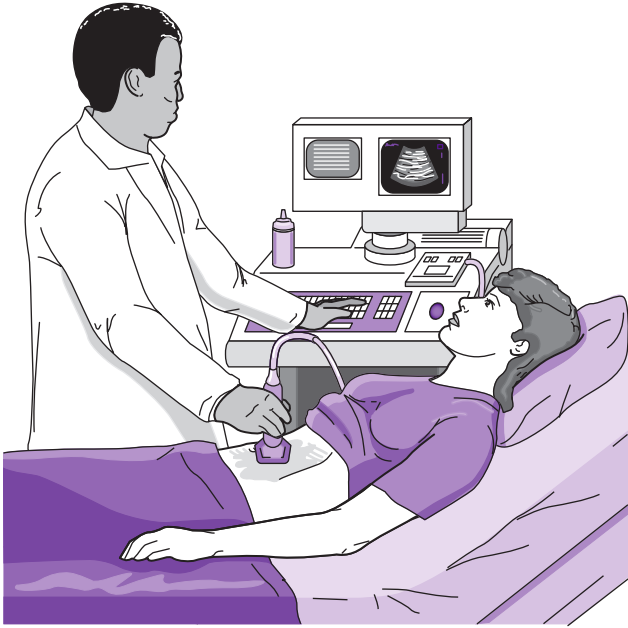

**FIG. 1** Ultrasound of the abdomen.

study requires no contrast material and has no associated radiation, it is especially useful in patients who are allergic to contrast, have diminished renal function, or are pregnant. Fasting may be preferred, but it is not mandatory.

### **Interfering factors**

- Barium blocks transmission of ultrasonic waves. For this reason, ultrasonography of the abdomen should be performed before any barium contrast studies.
- Large amounts of gas in the bowel distort visualization of abdominal organs because bowel gas reflects sound. Likewise, ultrasonic evaluation of the lungs yields poor results.
- Obesity may affect the results of the study because sound waves are altered by fatty tissue.
- The quality of the ultrasound image and the sufficiency of the study depend to a very large part on the abilities of the ultrasound technologist performing the study.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that fasting may or may not be required, depending on the organ to be examined. No fasting is required for ultrasonography of the abdominal aorta, kidney, liver, spleen, or pancreas. Fasting, however, is preferred for ultrasound of the gallbladder and bile ducts.

### During

- Note the following procedural steps:
  1. The patient is placed on the ultrasonography table in the prone or supine position, depending on the organ to be studied.
  2. A greasy conductive paste (coupling agent) is applied to the patient's skin. This paste is used to enhance sound wave transmission and reception.
  3. A transducer is placed over the skin.
  4. Pictures are taken of the reflections from the organs.
- The test is completed in about 1 hour, usually by an ultrasound technologist, and is interpreted by a radiologist.
- PT** Tell the patient that this procedure causes no discomfort.

### After

- Remove the coupling agent from the patient's skin.
- Note that if a biopsy is done, refer to biopsy of the specific organ (e.g., liver or kidney biopsy).

## Abnormal findings

### Kidney

Renal cysts  
 Renal tumor  
 Renal calculi  
 Hydronephrosis  
 Ureteral obstruction  
 Perirenal abscess  
 Glomerulonephritis  
 Pyelonephritis  
 Perirenal hematoma

### Gallbladder

Polyps  
 Tumor  
 Gallstone

## **4 abdominal ultrasound**

### **Liver**

Tumor

Abscess

Intrahepatic dilated bile ducts

### **Pancreas**

Tumor

Cysts

Pseudocysts

Abscess

Inflammation

### **Bile ducts**

Gallstone

Dilation

Stricture

Tumor

### **Abdominal aorta**

Aneurysm

### **Abdominal cavity**

Ascites

Abscess

### **notes**



## acetylcholine receptor antibody panel (AChR Ab, Anti-AChR antibody)

A

**Type of test** Blood

### Normal findings

ACh receptor (muscle) binding antibodies:  $\leq 0.02$  nmol/L

ACh receptor (muscle) modulating antibodies: 0%-20% (reported as % loss of AChR)

Striational (striated muscle) antibodies:  $< 1:60$

### Test explanation and related physiology

These antibodies may cause blocks in neuromuscular transmission by interfering with the binding of *acetylcholine* (ACh) to *ACh receptor* (AChR) sites on the muscle membrane, thereby preventing muscle contraction. This phenomenon characterizes myasthenia gravis (MG). Antibodies to AChR occur in more than 85% of patients with acquired MG. Lower levels are seen in patients with ocular MG only. The presence of these antibodies is virtually diagnostic of MG, but a negative test result does not exclude the disease. The measured titers do not correspond well with the severity of MG in different patients. In an individual patient, however, antibody levels are particularly useful in monitoring response to therapy. As the patient improves, antibody titers decrease. In adults with MG, there is at least a 20% occurrence of thymoma or other neoplasm. Neoplasms are an endogenous source of the antigens driving production of AChR autoantibodies.

Several AChR antibodies can be associated with MG. The AChR-*binding* antibody can activate complement and lead to loss of AChR. The AChR-*modulating* antibody causes receptor endocytosis, resulting in loss of AChR expression, which correlates most closely with clinical severity of disease. It is the most sensitive test. A positive modulating antibody test result may indicate subclinical MG, contraindicating the use of curare-like drugs during surgery. The AChR-*blocking* antibody may impair binding of ACh to the receptor, leading to poor muscle contraction. It is the least sensitive test (positive in only 61% of patients with MG).

*Antistriated muscle antibody (striated muscle antibody, IgG)* titers greater than or equal to 1:80 are suggestive of myasthenia. This antibody is detectable in 30% to 40% of anti-AChR-negative patients (particularly those with bulbar symptoms only). However, striated muscle antibody can be found in rheumatic fever, myocardial infarction, and a variety of postcardiotomy states.

## 6 acetylcholine receptor antibody panel

### Interfering factors

- False-positive results may occur in patients with amyotrophic lateral sclerosis who have been treated with cobra venom.
- False-positive results may be seen in patients with penicillamine-induced or Lambert-Eaton myasthenic syndrome.
- Patients with autoimmune liver disease may have elevated results.
- ▼ Drugs that may cause *increased* levels include muscle paralytic medicines (succinylcholine) and snake venom.
- ▼ Immunosuppressive drugs may suppress the formation of these antibodies in patients with subclinical MG.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased titer levels

Myasthenia gravis

Ocular myasthenia gravis

Thymoma

notes

## acid phosphatase (Prostatic acid phosphatase [PAP], Tartrate-resistant acid phosphatase [TRAP])

**Type of test** Blood

### Normal findings

Adult/elderly: 0.13-0.63 units/L (Roy, Brower, Hayden; 37° C)  
or 2.2-10.5 units/L (SI units)

Child: 8.6-12.6 units/mL (30° C)

Newborn: 10.4-16.4 units/mL (30° C)

### Test explanation and related physiology

Elevated levels are seen in patients with prostatic cancer that has metastasized beyond the capsule to other parts of the body, especially bone. The degree of elevation indicates the extent of disease. This test is rarely performed for this indication; there are better tumor markers (p. 194), including prostate-specific antigen (p. 743).

Because acid phosphatase is also found at high concentrations in seminal fluid, this test can be performed on vaginal secretions to investigate alleged rape. This is now the primary use of PAP testing. Acid phosphatase is a lysosomal enzyme. Therefore lysosomal storage diseases (e.g., Gaucher disease and Niemann-Pick disease) are associated with elevated levels.

### Interfering factors

- Alkaline and acid phosphatases are very similar enzymes that differ in the pH at which they are identified. Any condition associated with very high levels of alkaline phosphatase may falsely indicate high acid phosphatase levels.
- Falsely high levels of acid phosphatase may occur in males after a digital examination or after instrumentation of the prostate (e.g., cystoscopy) because of prostatic stimulation.
- ✚ Drugs that may cause *increased* levels of acid phosphatase include alglucerase, androgens (in females), and clofibrate.
- ✚ Drugs that may cause *decreased* levels include alcohol, fluorides, heparin, oxalates, and phosphates.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Avoid hemolysis. Red blood cells contain acid phosphatase.
- Note on the laboratory slip if the patient has had a prostatic examination or instrumentation of the prostate within the past 24 hours.

## 8 acid phosphatase

- Do *not* leave the specimen at room temperature for 1 hour or longer because the enzyme is heat and pH sensitive, and its activity will decrease.

### Abnormal findings

#### ▲ Increased levels

Prostatic carcinoma  
Benign prostatic hypertrophy  
Prostatitis  
Multiple myeloma  
Paget disease  
Hyperparathyroidism  
Metastasis to the bone  
Sickle cell crisis  
Thrombocytosis  
Lysosomal disorders (e.g., Gaucher disease)  
Renal diseases  
Liver diseases (e.g., cirrhosis)  
Rape

notes

**activated clotting time (ACT, Activated coagulation time)**

**Type of test** Blood

**Normal findings**

70-120 sec

Therapeutic range for anticoagulation: 150-600 sec

(Normal ranges and anticoagulation ranges vary according to type of laboratory procedure and particular therapy.)

**Possible critical values**

Depend on use for the test and clinical situation

**Test explanation and related physiology**

The ACT is primarily used to measure the anticoagulant effect of heparin or other direct thrombin inhibitors during cardiac angioplasty, hemodialysis, and cardiopulmonary bypass (CPB) surgery. This test measures the time for whole blood to clot after the addition of particulate activators. It is similar to the *activated partial thromboplastin time* (APTT, p. 681) in that it measures the ability of the *intrinsic* pathway to begin clot formation by activating factor XII (see [Figure 10](#), p. 251). By checking the blood clotting status with ACT, the response to unfractionated heparin therapy can be monitored.

Both the APTT and the ACT can be used to monitor heparin therapy for patients during CPB. However, the ACT has several advantages over the APTT. First, the ACT is more accurate than the APTT when high doses of heparin are used for anticoagulation. This makes it especially useful during clinical situations requiring high-dose heparin, such as during CPB, when high-dose anticoagulation is necessary at levels 10 times those used for venous thrombosis. The APTT is not measurable at these high doses. The accepted goal for the ACT is 400 to 480 seconds during CPB.

Second, the ACT is both less expensive and more easily performed, even at the bedside. This allows for immediate accessibility and decreased turnaround time. The capability to perform the ACT at the point of care makes the ACT particularly useful for patients requiring angioplasty, hemodialysis, and CPB.

A nomogram is often used as a guide to reach the desired level of anticoagulation. This nomogram is used in determining the dose of protamine to neutralize the heparin upon completion of these procedures. The ACT is used in determining when it is safe to remove the vascular access upon completion of these

## 10 activated clotting time

procedures. The benefits of the *modified* ACT test are that it requires a smaller-volume blood specimen; it can be automated; it can use standardized blood/reagent mixing; and it provides faster clotting time results than the conventional ACT. The modified ACT is now used more frequently.

### Interfering factors

- The ACT is affected by biologic variables, including hypothermia, hemodilution, and platelet number and function.
- Factors affecting the pharmacokinetics of heparin (e.g., kidney or liver disease) and heparin resistance can affect ACT measurements.
- A clotted specimen can increase ACT measurements.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: verify with laboratory
- Less than 1 mL of blood is collected and placed in a machine at the bedside. When a clot forms, the ACT value is displayed.
- If the patient is receiving a continuous heparin drip, the blood sample is obtained from the arm without the intravenous catheter.
- The bleeding time will be prolonged because of anticoagulation therapy.
- Assess the patient to detect possible bleeding. Check for blood in the urine and all other excretions, and assess the patient for bruises, petechiae, and low back pain.

### Abnormal findings

#### ▲ Increased levels

Heparin administration  
Clotting factor deficiencies  
Cirrhosis of the liver  
Lupus inhibitor  
Warfarin administration

#### ▼ Decreased levels

Thrombosis

notes

**adrenocorticotrophic hormone (ACTH, Corticotropin)**

**Type of test** Blood

**Normal findings**

Adult/elderly:

Female: 19 years and older: 6-58 pg/mL

Male: 19 years and older: 7-69 pg/mL

Children:

Male and female: 10-18 years: 6-55 pg/mL

Male and female: 1 week-9 years: 5-46 pg/mL

**Test explanation and related physiology**

The ACTH tests the anterior pituitary gland function and provides the greatest insight into the causes of either Cushing syndrome (overproduction of cortisol) or Addison disease (underproduction of cortisol). An elaborate feedback mechanism for cortisol exists to coordinate the function of the hypothalamus, pituitary gland, and adrenal glands. ACTH is an important part of this mechanism. Corticotropin-releasing hormone (CRH) is made in the hypothalamus. This stimulates ACTH production in the anterior pituitary gland. This, in turn, stimulates the adrenal cortex to produce cortisol. The rising levels of cortisol act as negative feedback and curtail further production of CRH and ACTH.

In a patient with Cushing syndrome, an elevated ACTH level can be caused by a pituitary or a nonpituitary (ectopic) ACTH-producing tumor, usually in the lung, pancreas, thymus, or ovary. ACTH levels over 200 pg/mL usually indicate ectopic ACTH production. If the ACTH level is lower than normal in a patient with Cushing syndrome, an adrenal adenoma or carcinoma is probably the cause of the hyperfunction.

In patients with Addison disease, an elevated ACTH level indicates primary adrenal gland failure, as in adrenal gland destruction caused by infarction, hemorrhage, or autoimmunity; surgical removal of the adrenal gland; congenital enzyme deficiency; or adrenal suppression after prolonged ingestion of exogenous steroids. If the ACTH level is lower than normal in a patient with adrenal insufficiency, hypopituitarism is most probably the cause of the hypofunction.

One must be aware that there is a diurnal variation of ACTH levels that corresponds to variation of cortisol levels. Levels in evening (8 PM to 10 PM) samples are usually one-half to two-thirds those of morning (4 AM to 8 AM) specimens. This diurnal

## 12 adrenocorticotrophic hormone

variation is lost when disease (especially neoplasm) affects the pituitary or adrenal glands. Likewise, stress can blunt or eliminate this normal diurnal variation.

### Interfering factors

- Stress (trauma, pyrogens, or hypoglycemia) and pregnancy can increase levels.
- Recently administered radioisotope scans can affect levels.
- ▼ Drugs that may cause *increased* ACTH levels include aminoglutethimide, amphetamines, estrogens, ethanol, insulin, metyrapone, spironolactone, and vasopressin.
- ▼ Corticosteroids may *decrease* ACTH levels.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: green
- Evaluate the patient for stress factors that could invalidate the test results.
- Evaluate the patient for sleep pattern abnormalities. With a normal sleep pattern, the ACTH level is the highest between 4 AM and 8 AM and the lowest around 9 PM.
- Chill the blood tube to prevent enzymatic degradation of ACTH.
- Place the specimen in ice water and send it to the chemistry laboratory immediately. ACTH is a very unstable peptide in plasma and should be stored at  $-20^{\circ}\text{C}$  to prevent artificially low values.

### Abnormal findings

#### ▲ Increased levels

Addison disease (primary adrenal insufficiency)  
Cushing syndrome (pituitary-dependent adrenal hyperplasia)  
Ectopic ACTH syndrome  
Stress  
Adrenogenital syndrome (congenital adrenal hyperplasia)

#### ▼ Decreased levels

Secondary adrenal insufficiency (pituitary insufficiency)  
Cushing syndrome  
Hypopituitarism  
Adrenal adenoma or carcinoma  
Steroid administration



## adrenocorticotropic hormone stimulation test with cosyntropin (ACTH stimulation test, Cortisol stimulation test)

**Type of test** Blood

### Normal findings

Rapid test: cortisol levels increase  $> 7$  mcg/dL higher than baseline

24-hour test: cortisol levels  $> 40$  mcg/dL

3-day test: cortisol levels  $> 40$  mcg/dL

### Test explanation and related physiology

This test is performed in patients found to have an adrenal insufficiency. An increase in plasma cortisol levels after the infusion of an ACTH-like drug indicates that the adrenal gland is normal and is capable of functioning if stimulated. In that case, the cause of the adrenal insufficiency would lie within the pituitary gland (hypopituitarism, which is called secondary adrenal insufficiency). If little or no rise in cortisol levels occurs after the administration of the ACTH-like drug, the adrenal gland is the source of the problem and cannot secrete cortisol. This is called primary adrenal insufficiency (Addison disease), which may be caused by adrenal hemorrhage, infarction, autoimmunity, metastatic tumor, surgical removal of the adrenal glands, or congenital adrenal enzyme deficiency.

This test can also be used in the evaluation of patients with Cushing syndrome. Patients with Cushing syndrome caused by bilateral adrenal hyperplasia have an exaggerated cortisol elevation in response to the administration of the ACTH-like drug. Those experiencing Cushing syndrome as a result of hyperfunctioning adrenal tumors (which are usually autonomous and relatively insensitive to ACTH) have little or no increase in cortisol levels over baseline values.

Cosyntropin (Cortrosyn) is a synthetic subunit of ACTH that has the same corticosteroid-stimulating effect as endogenous ACTH in healthy persons. During this test, cosyntropin is administered to the patient, and the ability of the adrenal gland to respond is measured by plasma cortisol levels.

The *rapid stimulation test* is only a screening test. A normal response excludes adrenal insufficiency. An abnormal response, however, requires a 24-hour to 3-day prolonged ACTH stimulation test to differentiate primary insufficiency from secondary insufficiency. It should be noted that the adrenal gland can also be stimulated by insulin-induced hypoglycemia as a stressing agent. When insulin is the stimulant, cortisol and glucose levels are measured.

## 14 ACTH stimulation test with cosyntropin

### Interfering factors

- ▼ Drugs that may cause artificially *increased* cortisol levels include corticosteroids, estrogens, and spironolactone.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red

### Rapid test

- Obtain a baseline plasma cortisol level. This should be done 30 minutes before cosyntropin (ACTH-like drug) administration.
- Administer an IV injection of cosyntropin over a 2-minute period as prescribed.
- Measure plasma cortisol levels 30 and 60 minutes after drug administration.

### 24-hour test

- Obtain a baseline plasma cortisol level.
- Start an IV infusion of synthetic cosyntropin.
- Administer the solution as prescribed for 24 hours.
- After 24 hours, obtain another plasma cortisol level.

### 3-day test

- Obtain a baseline plasma cortisol level.
- Administer the prescribed dose of cosyntropin IV over an 8-hour period for 2 to 3 consecutive days.
- Measure plasma cortisol levels at 12, 24, 36, 48, 60, and 72 hours after the start of the test.

### Abnormal findings

#### In adrenal insufficiency

*Increase higher than normal response (secondary adrenal insufficiency)*

Hypopituitarism

Exogenous steroid ingestion

Endogenous steroid production from a nonendocrine tumor

*Normal or lower than normal response (primary adrenal insufficiency)*

Addison disease

Adrenal infarction/hemorrhage

Metastatic tumor to the adrenal gland

Congenital enzyme adrenal insufficiency

Surgical removal of the adrenal gland

**In Cushing syndrome**

*Increase higher than normal response*

Bilateral adrenal hyperplasia

*Normal or lower than normal response*

Adrenal adenoma

Adrenal carcinoma

ACTH-producing nonadrenal tumor

Chronic steroid use

**notes**

**adrenocorticotrophic hormone stimulation test with metyrapone** (ACTH stimulation test with metyrapone, Metyrapone test)

**Type of test** Blood; urine (24-hour)

**Normal findings**

**Blood**

11-deoxycortisol increased to  $> 7$  mcg/dL and cortisol  $< 10$  mcg/dL

**Urine (24-hour)**

Baseline excretion of urinary *17-hydroxycorticosteroid* (17-OCHS) more than doubled

**Test explanation and related physiology**

Metyrapone is a potent blocker of an enzyme involved in cortisol production. Therefore cortisol production is reduced. When this drug is given, the resulting fall in cortisol production should stimulate pituitary secretion of ACTH by way of a negative feedback mechanism. Cortisol precursors (*11-deoxycortisol* and *17-OCHS*) can be detected in the urine or blood. This test is similar to the ACTH stimulation test with *cosyntropin* (see p. 13).

In patients with adrenal hyperplasia caused by pituitary overproduction of ACTH, the cortisol precursors are greatly increased, more than expected in normal patients. This is because the normal adrenal-pituitary response mechanism is still intact. No response to metyrapone occurs in patients with Cushing syndrome resulting from adrenal adenoma or carcinoma because the tumors are autonomous and therefore insensitive to changes in ACTH secretion.

This test is also used to evaluate the pituitary reserve capacity to produce ACTH. It can document that adrenal insufficiency exists as a result of pituitary disease (secondary adrenal insufficiency) rather than primary adrenal pathology.

**Contraindications**

- Patients with possible adrenal insufficiency
- Patients taking glucocorticosteroids

**Potential complications**

- Addison disease and Addisonian crisis because metyrapone inhibits cortisol production
- Dizziness, sedation, allergic reaction, and bone marrow suppression

## Interfering factors

- Recent administration of radioisotopes can affect results.
- ✚ Chlorpromazine interferes with the response to metyrapone and should not be administered during the testing.

## Procedure and patient care

- Obtain a baseline cortisol level (see p. 290) for the blood test.
- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or green
- See inside front cover for Routine Urine Testing.
- Obtain a baseline 24-hour urine specimen for the 17-OCHS level (see p. 526) for the urine test.

## During

### *Blood*

- Administer a prescribed dose of metyrapone at 11 PM the night before the blood sample is to be collected. Collect a venous blood sample in a red-top tube in the morning.

### *Urine*

- Obtain a 24-hour urine specimen for the 17-OCHS level as a baseline. Then collect a 24-hour urine specimen for the 17-OCHS level during and again 1 day after the oral administration of a dose of metyrapone, which is given in 5 doses every 4 hours within 24 hours.
- Metyrapone should be administered with a glass of milk to diminish any GI side effects.

## After

- Assess the patient for impending signs of Addisonian crisis (muscle weakness, mental and emotional changes, anorexia, nausea, vomiting, hypotension, hyperkalemia, or vascular collapse).
- Note that Addisonian crisis is a medical emergency that must be treated vigorously with replenishing steroids, reversing shock, and restoring circulation.

## Abnormal findings

### Increased cortisol precursors

Adrenal hyperplasia

### No change in cortisol precursors

Adrenal tumor

Ectopic ACTH syndrome

Secondary adrenal insufficiency

## notes

**adrenal steroid precursors** (Androstenediones [AD], Dehydroepiandrosterone [DHEA], Dehydroepiandrosterone sulfate [DHEA S], 11-Deoxycortisol, 17-Hydroxyprogesterone, 17-Hydroxypregnenolone, Pregnenolone)

**Type of test** Blood

**Normal findings**

		Female	Male
AD	Tanner stage I	0.05-0.51 ng/mL	0.04-0.32 ng/mL
	Tanner stage II	0.15-1.37 ng/mL	0.08-0.48 ng/mL
	Tanner stage III	0.37-2.24 ng/mL	0.14-0.87 ng/mL
	Tanner stage IV-V	0.35-2.05 ng/mL	0.27-1.07 ng/mL
DHEA	Tanner stage I	0.14-2.76 ng/mL	0.11-2.37 ng/mL
	Tanner stage II	0.83-4.87 ng/mL	0.37-3.66 ng/mL
	Tanner stage III	1.08-7.56 ng/mL	0.75-5.24 ng/mL
	Tanner stage IV-V	1.24-7.88 ng/mL	1.22-6.73 ng/mL
DHEA S	Tanner stage I	7-209 mcg/dL	7-126 mcg/dL
	Tanner stage II	28-260 mcg/dL	13-241 mcg/dL
	Tanner stage III	39-390 mcg/dL	32-446 mcg/dL
	Tanner stage IV-V	81-488 mcg/dL	65-371 mcg/dL

**Test explanation and related physiology**

Androstenediones (ADs, DHEA, and the sulfuric ester, DHEA S) are precursors of testosterone and estrone and are made in the gonads and the adrenal gland. 11-Deoxycortisol, 17-hydroxyprogesterone, 17-hydroxypregnenolone, and pregnenolone are precursors of cortisol. ACTH stimulates their adrenal secretion. Children with congenital adrenal hyperplasia (CAH) have genetic mutations that cause deficiencies in the enzymes involved in the synthesis of cortisol, testosterone, aldosterone, and estrone. When defects in enzyme synthesis occur along the path of hormone synthesis, the previously listed precursors exist in levels that exceed normal through the increased stimulation of ACTH.

The symptoms of the disorder depend on which steroids are overproduced and which are deficient. As a result, CAH may present with various symptoms, including virilization of the affected female infant, signs of androgen excess in males and females, signs of sex hormone deficiency in males and females, salt-wasting crisis secondary to cortisol and aldosterone deficiency, or hormonal hypertension due to increased mineralocorticoids. A milder, nonclassic form of CAH is characterized by premature puberty, acne, hirsutism, menstrual irregularity, and infertility.

These same precursors can occur in adults due to adrenal or gonadal tumors. Patients with polycystic ovary syndrome (Stein-Leventhal syndrome) have particularly elevated levels of ADs. DHEA S levels are particularly high in patients with adrenal carcinoma.

In patients suspected of CAH, testing for a panel of steroids involved in the cortisol biosynthesis pathway may be performed to establish the specific enzyme deficiency. In most cases, basal concentrations within the normal reference interval rule out CAH. The ratio of the precursor to the final pathway product (with and without ACTH stimulation) may be used to diagnose which enzyme is deficient.

### Interfering factors

- A radioactive scan performed 1 week before the test may invalidate the test results if radioimmunoassay is performed.
- ✚ Drugs that may *increase* levels of ADs include clomiphene, corticotropin, and metyrapone.
- ✚ Drugs that may *decrease* levels of ADs include steroids.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: serum separator or red
- PT Tell the female patient that the specimen should be collected 1 week before or after the menstrual period.
- Indicate the date of the last menstrual period (if applicable) on the laboratory form.

### Abnormal findings

#### ▲ Increased levels

Adrenal tumor  
 Congenital adrenal hyperplasia  
 Ectopic ACTH-producing tumors  
 Cushing syndrome (some cases)  
 Stein-Leventhal syndrome  
 Ovarian sex cord tumor

#### ▼ Decreased levels

Gonadal failure  
 Primary or secondary adrenal insufficiency

**age-related macular degeneration risk analysis** (ARMD risk analysis, Y402H, and A69S)**Type of test** Blood**Normal findings**

No mutation noted

**Test explanation and related physiology**

Age-related macular degeneration (ARMD) is recognized as a leading cause of blindness in the United States. Blurred or distorted vision and difficulty adjusting to dim light are common symptoms. ARMD, both wet and dry types, is considered a multifactorial disorder because it is thought to develop because of an interplay between environmental (smoking) and genetic (gender, ethnicity) risk and protective (antioxidants) factors. At least two genetic variants (Y402H and A69S) have been found to be associated with an increased risk for ARMD. The Y402H and the A69S genetic variants are common polymorphisms in ARMD. An individual with two copies of the Y402H variant in the gene *CFH* and two copies of the A69S variant in the gene *LOC387715* has an approximately 60-fold increased risk for ARMD. This is significant, given how common ARMD is in the general population.

This information can be clinically useful when making medical management decisions (e.g., the use of inflammatory markers) and emphasizing to patients the benefits of smoking cessation and dietary modification. In some cases, genotype information may also assist with clinical diagnosis.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
  - Fasting: no
  - Blood tube commonly used: lavender or yellow
- PT** Tell the patient that results may not be available for a few weeks.

**Abnormal findings**

Increased risk of ARMD

**notes**



**alanine aminotransferase** (ALT, formerly Serum glutamic-pyruvic transaminase [SGPT])

**Type of test** Blood

**Normal findings**

Adult/child: 4-36 units/L at 37° C, or 4-36 units/L (SI units)

Elderly: may be slightly higher than adult

Infant: may be twice as high as adult

**Test explanation and related physiology**

ALT is found predominantly in the liver; lesser quantities are found in the kidneys, heart, and skeletal muscles. Injury or disease affecting the liver parenchyma causes a release of this hepatocellular enzyme into the bloodstream, thus elevating serum ALT levels. Generally, most ALT elevations are caused by liver disease. Therefore this enzyme is not only sensitive but also very specific in indicating hepatocellular disease. In hepatocellular disease other than viral hepatitis, the ALT/AST ratio (*De Ritis ratio*) is less than 1. In viral hepatitis, the ratio is greater than 1. This is helpful in the diagnosis of viral hepatitis.

**Interfering factors**

- Previous IM injections may cause elevated levels.
- ✱ Drugs that may cause *increased* ALT levels include acetaminophen, allopurinol, aminosalicic acid (PAS), ampicillin, azathioprine, carbamazepine, cephalosporins, chlordiazepoxide, chlorpropamide, clofibrate, cloxacillin, codeine, dicumarol, indomethacin, isoniazid (INH), methotrexate, methyldopa, nafcillin, nalidixic acid, nitrofurantoin, oral contraceptives, oxacillin, phenothiazines, phenylbutazone, phenytoin, procainamide, propoxyphene, propranolol, quinidine, salicylates, tetracyclines, and verapamil.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Patients with liver dysfunction often have prolonged clotting times.

## 22 alanine aminotransferase

### Abnormal findings

#### ▲ Increased levels

- Hepatitis
- Hepatic necrosis
- Hepatic ischemia
- Cirrhosis
- Cholestasis
- Hepatic tumor
- Hepatotoxic drugs
- Obstructive jaundice
- Severe burns
- Trauma to striated muscle
- Myositis
- Pancreatitis
- Myocardial infarction
- Infectious mononucleosis
- Shock

#### notes

## aldolase

**Type of test** Blood

### Normal findings

Adult: 3-8.2 Sibley-Lehninger units/dL or 22-59 mU/L at 37° C (SI units)

Child: approximately two times the adult values

Newborn: approximately four times the adult values

### Test explanation and related physiology

Serum aldolase is very similar to the enzymes aspartate aminotransferase AST (SGOT) (see p. 125) and CPK (see p. 297). Aldolase is an enzyme used in glycolysis (breakdown of glucose). As with AST and creatine phosphokinase, aldolase exists throughout the body in most tissues. This test is most useful for indicating muscular or hepatic cellular injury or disease. The serum aldolase level is very high in patients with muscular dystrophies, dermatomyositis, and polymyositis. Levels also are increased in patients with gangrenous processes, muscular trauma, and muscular infectious diseases (e.g., trichinosis). Elevated levels are also noted in chronic hepatitis, obstructive jaundice, and cirrhosis.

Neurologic diseases causing weakness can be differentiated from muscular causes of weakness with this test. Normal values are seen in patients with such neurologic diseases as poliomyelitis, myasthenia gravis, and multiple sclerosis. Elevated aldolase levels are seen in the primary muscular disorders.

### Interfering factors

- Previous IM injections may cause elevated levels.
- Strenuous exercise can cause a transient spike in aldolase.
- ✚ Drugs that may cause *increased* aldolase levels include hepatotoxic agents.
- ✚ Drugs that may cause *decreased* aldolase levels include phenothiazines.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: verify with laboratory
- Blood tube commonly used: red

## 24 aldolase

### Abnormal findings

#### ▲ Increased levels

- Hepatocellular diseases (e.g., hepatitis)
- Muscular diseases (e.g., muscular dystrophy, dermatomyositis, and polymyositis)
- Muscular trauma (e.g., severe crush injuries)
- Muscular infections (e.g., trichinosis)
- Gangrenous processes (e.g., gangrene of the bowel)
- Myocardial infarction

#### ▼ Decreased levels

- Late muscular dystrophy
- Hereditary fructose intolerance
- Muscle-wasting disease

### notes

## aldosterone

A

**Type of test** Blood; urine (24-hour)

### Normal findings

#### Blood

Supine: 3-10 ng/dL or 0.08-0.30 nmol/L (SI units)

Upright:

Female: 5-30 ng/dL or 0.14-0.80 nmol/L (SI units)

Male: 6-22 ng/dL or 0.17-0.61 nmol/L (SI units)

Child/adolescent:

Newborn: 5-60 ng/dL

1 week-1 year: 1-160 ng/dL; 5-7 years: 5-50 ng/dL

1-3 years: 5-60 ng/dL; 7-11 years: 5-70 ng/dL

3-5 years: 5-80 ng/dL; 11-15 years: 5-50 ng/dL

#### Urine (24-hour)

2-26 mcg/24 hr or 6-72 nmol/24 hr (SI units)

### Test explanation and related physiology

This test is used to diagnose hyperaldosteronism. Production of aldosterone, a hormone produced by the adrenal cortex, is regulated primarily by the renin-angiotensin system. Secondarily, aldosterone is stimulated by ACTH, low serum sodium levels, and high serum potassium levels. Aldosterone in turn stimulates the renal tubules to absorb sodium (water follows) and to secrete potassium into the urine. In this way, aldosterone regulates serum sodium and potassium levels. Because water follows sodium transport, aldosterone also partially regulates water absorption (and plasma volume).

Increased aldosterone levels are associated with primary aldosteronism, in which a tumor (usually an adenoma) of the adrenal cortex (Conn syndrome) or bilateral adrenal nodular hyperplasia causes increased production of aldosterone. Patients with primary aldosteronism characteristically have hypertension, weakness, polyuria, and hypokalemia.

Increased aldosterone levels also occur with secondary aldosteronism caused by nonadrenal conditions. These include:

- Renal vascular stenosis or occlusion
- Hyponatremia (from diuretic or laxative abuse) or low salt intake
- Hypovolemia
- Pregnancy or use of estrogens
- Malignant hypertension
- Potassium loading

## 26 aldosterone

- Edematous states (e.g., congestive heart failure, cirrhosis, nephrotic syndrome)

The aldosterone assay can be done on a 24-hour urine specimen or a plasma blood sample. The advantage of the 24-hour urine sample is that short-term fluctuations are eliminated. Plasma values are more convenient to sample, but they are affected by the short-term fluctuations.

Primary aldosteronism can be diagnosed by demonstrating very little to no rise in serum renin levels after an *aldosterone stimulation test* (using salt restriction as the stimulant). This is because aldosterone is already maximally secreted by the pathologic adrenal gland. Failure to suppress aldosterone with saline infusion (1.5-2 L of NSS infused between 8 PM and 10 AM, called an *aldosterone suppression test*) is further evidence of primary aldosteronism. Aldosterone can also be measured in blood obtained from adrenal venous sampling.

### Interfering factors

- Strenuous exercise and stress can stimulate adrenocortical secretions and increase aldosterone levels.
- Excessive licorice ingestion can cause decreased levels because it produces an aldosterone-like effect.
- Values are influenced by posture, position, diet, diurnal variation, and pregnancy.
- If the test is performed using radioimmunoassay, recently administered radioactive medications will affect test results.
- ✚ Drugs that may cause *increased* levels include diazoxide, diuretics, hydralazine, laxatives, nitroprusside, potassium, and spironolactone.
- ✚ Drugs that may cause *decreased* levels include angiotensin-converting inhibitors (e.g., captopril), fludrocortisone, and propranolol, as well as licorice.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: serum separator
- PT** Note that the patient is asked to be in the upright position (at least sitting) for at least 2 hours before the blood is drawn.
- PT** Explain the procedure for collecting a 24-hour urine sample if urinary aldosterone is ordered. (See inside front cover for Routine Urine Testing.)
- PT** Give the patient verbal and written instructions regarding dietary and medication restrictions.

- PT** Instruct the patient to maintain a normal sodium diet (~3 g/day) for at least 2 weeks before the blood or urine collection.
- PT** Have the patient ask the physician whether drugs that alter sodium, potassium, and fluid balance (e.g., diuretics, antihypertensives, steroids, oral contraceptives) should be withheld. Test results will be more accurate if these are suspended at least 2 weeks before either the blood or the urine test.
- PT** Inform the patient that renin inhibitors (e.g., propranolol) should not be taken 1 week before the test.
- PT** Tell the patient to avoid licorice for at least 2 weeks before the test because of its aldosterone-like effect.

### During

- Occasionally, for hospitalized patients, draw the sample with the patient in the supine position *before* he or she rises.
- Obtain the specimen in the morning.
- Note that sometimes a second specimen (upright sample) is collected 4 hours later, after the patient has been up and moving.

### After

- Indicate on the laboratory slip if the patient was supine or standing during the venipuncture.
- Handle the blood specimen gently. Rough handling may cause hemolysis and alter the test results.
- Transport the specimen on ice to the laboratory.

### Abnormal findings

#### ▲ Increased levels

##### *Primary aldosteronism*

Aldosterone-producing adrenal adenoma (Conn syndrome)  
 Adrenal cortical nodular hyperplasia  
 Bartter syndrome

##### *Secondary aldosteronism*

Hyponatremia  
 Hyperkalemia  
 Diuretic ingestion resulting in hypovolemia and hyponatremia  
 Laxative abuse

#### ▼ Decreased levels

Aldosterone deficiency  
 Renin deficiency  
 Steroid therapy  
 Addison disease  
 Patients on a high sodium diet  
 Hyponatremia  
 Hypokalemia  
 Toxemia of pregnancy  
 Antihypertensive therapy

## 28 aldosterone

### ▲ Increased levels

*Secondary aldosteronism*

(*cont.*)

Stress

Malignant hypertension

Generalized edema

Renal arterial stenosis

Pregnancy

Oral contraceptives

Hypovolemia or

hemorrhage

Cushing syndrome

notes



## alkaline phosphatase (ALP)

**Type of test** Blood

### Normal findings

Adult: 30-120 units/L or 0.5-2.0  $\mu$ Kat/L

Elderly: slightly higher than adults

Child/adolescent:

< 2 years: 85-235 units/L

2-8 years: 65-210 units/L

9-15 years: 60-300 units/L

16-21 years: 30-200 units/L

### Test explanation and related physiology

Although ALP is found in many tissues, the highest concentrations are found in the liver, biliary tract epithelium, and bone. Detection of this enzyme is important for determining liver and bone disorders. Enzyme levels of ALP are greatly increased in both extrahepatic and intrahepatic obstructive biliary disease and cirrhosis. Reports have indicated that the most sensitive test to indicate metastatic tumor to the liver is ALP.

Bone is the most frequent extrahepatic source of ALP; new bone growth is associated with elevated ALP levels, which explains why ALP levels are high in adolescents. Pathologic new bone growth occurs with osteoblastic metastatic (e.g., breast, prostate) tumors. Paget disease, healing fractures, rheumatoid arthritis, hyperparathyroidism, and normal-growing bones are sources of elevated ALP levels as well.

*Isoenzymes* of ALP are sometimes used to distinguish between liver and bone diseases. The detection of isoenzymes can help differentiate the source of the pathology associated with the elevated total ALP. ALP<sub>1</sub> is from the liver. ALP<sub>2</sub> is from the bone.

### Interfering factors

- Recent ingestion of a meal can increase ALP levels.
- ✚ Drugs that may cause *elevated* ALP levels include albumin made from placental tissue, allopurinol, antibiotics, azathioprine, colchicine, fluorides, indomethacin, isoniazid (INH), methotrexate, methyl dopa, nicotinic acid, phenothiazine, probenecid, tetracyclines, and verapamil.
- ✚ Drugs that may cause *decreased* levels include arsenicals, cyanides, fluorides, nitrofurantoin, oxalates, and zinc salts.

## 30 alkaline phosphatase

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Note that overnight fasting may be required for isoenzymes.
- Patients with liver dysfunction often have prolonged clotting times.

### Abnormal findings

#### ▲ Increased levels

Cirrhosis  
Intrahepatic or extrahepatic  
biliary obstruction  
Primary or metastatic liver  
tumor  
Intestinal ischemia or  
infarction  
Metastatic tumor to the  
bone  
Healing fracture  
Hyperparathyroidism  
Paget disease of bone  
Rheumatoid arthritis  
Sarcoidosis  
Osteomalacia  
Rickets

#### ▼ Decreased levels

Hypothyroidism  
Malnutrition  
Milk-alkali syndrome  
Pernicious anemia  
Hypophosphatemia  
Scurvy (vitamin C  
deficiency)  
Celiac disease  
Excess vitamin B  
ingestion  
Hypophosphatasia

### notes

**allergy blood testing** (IgE antibody test, Radioallergosorbent test [RAST])

**Type of test** Blood

**Normal findings**

**Total IgE serum**

Adult: 0-100 international units/mL

Child:

0-23 months: 0-13 international units/mL

2-5 years: 0-56 international units/mL

6-10 years: 0-85 international units/mL

**Test explanation and related physiology**

Measurement of serum IgE is an effective method to diagnose allergy and specifically identify the allergen (the substance to which the person is allergic). Serum IgE levels increase when allergic individuals are exposed to the allergen. Various classes of allergens can initiate the allergic response.

Although skin testing (see p. 33) can also identify a specific allergen, measurement of serum levels of IgE is helpful when a skin test result is questionable, when the allergen is not available in a form for dermal injection, or when the allergen may incite an anaphylactic reaction if injected. IgE is particularly helpful in cases in which skin testing is difficult (e.g., in infants or in patients with dermatographism or widespread dermatitis), and it is not always necessary to remove the patient from antihistamines, ACE inhibitors, antidepressants, or beta-blockers. The decision as to which method to use to diagnose an allergy and to identify the allergen depends on the elapsed time between exposure to an allergen and testing, class of allergen, the age of the patient, the possibility of anaphylaxis, and the affected target organ (e.g., skin, lungs, or intestine). In general, allergy skin testing is the preferred method in comparison with various *in vitro* tests for assessing the presence of specific IgE antibodies because it is more sensitive and specific, simpler to use, and less expensive.

IgE levels, similar to provocative skin testing, are used not only to diagnose allergy but also to identify the allergen so that an immunotherapeutic regimen can be developed. Increased levels of total IgE can be diagnostic of allergic disease in general. Specific IgE blood allergy testing, however, is an *in vitro* test for specific IgE directed to a specific allergen. Since the development of liquid allergen preparations, the use of *in vitro* blood allergy testing has increased considerably. It is more accurate and safer

## 32 allergy blood testing

than skin testing. Allergy testing of IgG antibodies can also be performed and may provide a more accurate correlation between allergen and allergic symptoms. Similar to IgE antibody testing, IgG antibody testing is often performed in “panels.” For example, there are meat panels that might include IgE or IgG testing for chicken, duck, goose, and turkey. Testing a fruit panel might include IgE or IgG antibody testing for apples, bananas, peaches, and pears. Testing in panels diminishes the cost of testing. Specific allergen antibody testing can follow panel testing.

### Contraindications

- Patients with multiple allergies; no information will be obtained regarding identification of the specific allergen.

### Interfering factors

- Concurrent diseases associated with elevated IgG levels will cause false-negative results.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: serum separator
- PT Inform the patient that the suspected allergen will be mixed with the patient’s blood specimen in the laboratory. The patient will not experience any allergic reaction by this method of testing.
- Determine whether the patient has recently been treated with a corticosteroid for allergies.

### Abnormal findings

Allergy-related diseases

Asthma

Dermatitis

Food allergy

Drug allergy

Occupational allergy

Allergic rhinitis

Angioedema

notes

## allergy skin testing

**Type of test** Skin

### Normal findings

< 3 mm wheal diameter

< 10 mm flare diameter

### Test explanation and related physiology

When properly performed, skin testing is the most convenient and least expensive test for detecting allergic reactions. Skin testing provides useful confirmatory evidence when a diagnosis of allergy is suspected on clinical grounds. The simplicity, rapidity, low costs, sensitivity, and specificity explain the crucial position skin testing has in allergy testing.

In an allergic patient, immediate wheal (swelling) and flare (redness) reactions follow injection of the specific allergen (that substance to which the person is allergic). This reaction is initiated by IgE and is mediated primarily by histamine secreted from mast cells. This usually occurs in about 5 minutes and peaks at 30 minutes. In some patients, a *late-phase reaction* occurs, which is highlighted by antibody and cellular infiltration into the area. This usually occurs within 1 to 2 hours.

There are three commonly accepted methods of injecting the allergen into the skin. The first method is called the *prick-puncture test* or *scratch test*. In this method, the allergen is injected into the epidermis. Life-threatening anaphylaxis reactions have not been reported with this method. The second method is called the *intradermal test*. Here the allergen is injected into the dermis (creating a skin wheal). Large local reactions and anaphylaxis have been reported with this latter method. For these two tests, the allergen placement part of the test takes about 5 to 10 minutes. The third method is called the *patch test*. This takes much longer because the patient must wear the patch for 48 hours to see if there is a delayed allergic reaction. With this method, needles are not used. Instead, an allergen is applied to a patch that is placed on the skin. It is usually done to detect whether a particular substance (e.g., latex, medications, fragrances, preservatives, hair dyes, metals, resins) is causing an allergic skin irritation, such as contact dermatitis.

Patients with dermographism (nonallergic response of redness and swelling of the skin at the site of any stimulation) develop a skin wheal with any skin irritation, even if nonallergic. In these patients, a false-positive reaction can occur with skin testing.

## 34 allergy skin testing

To eliminate these sorts of false positives, a “negative control” substance consisting of just the diluent without an allergen is injected at the same time as the other skin tests are performed. Patients who are immunosuppressed because of concurrent disease or medicines may have a blunted skin reaction even in the face of allergy. This would cause false-negative results. To avoid false negatives, a “positive control” substance consisting of a histamine analog is also injected into the forearm at the time of skin testing. This will cause a wheal and flare response even in the nonallergic patient unless the patient is immunosuppressed.

For inhalant allergens, skin tests are extremely accurate. However, for food allergies, latex allergies, drug sensitivity, and occupational allergies, skin tests are less reliable.

### Contraindications

- Patients with a history of prior anaphylaxis

### Potential complications

- Anaphylaxis

### Interfering factors

- False-positive results may occur with dermographism.
- False-positive results may occur if the patient has a reaction to the diluent used to preserve the extract.
- False-negative results may be caused by poor-quality allergen extracts, diseases that attenuate the immune response, or improper technique.
- Infants and the elderly may have decreased skin reactivity.
- ✶ Drugs that may *decrease* the immune response of skin testing include ACE inhibitors, beta-blockers, corticosteroids, nifedipine, and theophylline.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Observe skin testing precautions:
  1. Be sure that a physician is immediately available.
  2. Evaluate the patient for dermographism.
  3. Have medications and equipment available to handle anaphylaxis.
  4. Proceed with caution in patients with current allergic symptoms.
  5. Render great detail to the injection technique.
  6. Avoid bleeding due to injection.
  7. Avoid spreading of allergen solutions during the test.
  8. Record the skin reaction at the proper time.

- Obtain a history to evaluate the risk of anaphylaxis.
- Identify any immunosuppressive medications the patient may be taking.
- Evaluate the patient for dermographism by rubbing the skin with a pencil eraser and looking for a wheal at the site of irritation.
- Draw up 0.05 mL of 1:1000 aqueous epinephrine into a syringe before testing in the event of an exaggerated allergic reaction.
- A negative prick-puncture test should be performed before an intradermal test.

## During

### *Prick-puncture method*

- A drop of the allergen solution is placed onto the volar surface of the forearm or back.
- A 25-gauge needle is passed through the droplet and inserted into the epidermal space at angle with the bevel facing up.
- The skin is lifted up and the fluid is allowed to seep in. Excess fluid is wiped off after about a minute.

### *Intradermal method*

- With a 25-gauge needle, the allergen solution is injected into the dermis by creating a skin wheal. In this method, the bevel of the needle faces downward. A volume of between 0.01 and 0.05 mL is injected.
- In general, the allergen solution is diluted 100- to 1000-fold before injection.

### *Patch method*

- Clean the skin area (usually the back or arm).
- Apply the patches to the skin (as many as 20-30 can be applied).
- Instruct the patient to wear the patches for 48 hours. Tell the patient to avoid bathing or activities that cause heavy sweating.
- Tell the patient the patches will be removed at the doctor's office. Irritated skin at a patch site may indicate an allergy.

## After

- Evaluate the patient for an exaggerated allergic response.
- In the event of a systemic reaction, a tourniquet should be placed above the testing site, and epinephrine should be administered subcutaneously.
- With a pen, circle the area of testing and mark the allergen used.
- Read the skin test at the appropriate time.

## 36 allergy skin testing

- Skin tests are read when the reaction is mature, after about 15 to 20 minutes. Both the largest and smallest diameters of the wheal are determined. The measurements are averaged.
- The flare is measured in the same manner.
- Observe the patient for 20 to 30 minutes before discharge.

### Abnormal findings

Allergy-related diseases

Asthma

Dermatitis

Food allergy

Drug allergy

Occupational allergy

Allergic rhinitis

Angioedema

**notes**



**alpha<sub>1</sub>-antitrypsin** (A<sub>1</sub>,AT, AAT, Alpha<sub>1</sub>-antitrypsin phenotyping)

**Type of test** Blood

**Normal findings**

85-213 mg/dL or 0.85-2.13 g/L (SI units)

Negative for S or Z phenotype

**Test explanation and related physiology**

Serum alpha<sub>1</sub>-antitrypsin (AAT) levels should be obtained when an individual has a family history of early emphysema. AAT mutated genes can be associated with reduced levels of this enzyme. Patients with mutated AAT genes develop severe, disabling lower lung emphysema in the third and fourth decades of life. A similar deficiency is seen in children with cirrhosis and other liver diseases. The severity of these diseases is determined by the ATT phenotype.

Routine serum protein electrophoresis (see p. 746) is a good screening test for AAT deficiency because AAT accounts for about 90% of the protein in the alpha<sub>1</sub>-globulin region. Quantitative serum levels indicate deficiency. If levels are low, ATT phenotyping can be performed. The three major alleles determined by phenotyping are: M (full functioning, normal allele), S (associated with reduced levels of protein), and Z (disease-causing mutation associated with liver disease and premature emphysema). The S and Z alleles account for the majority of the abnormal alleles detected in affected patients.

*Inherited* AAT deficiency is associated with symptoms earlier in life than acquired AAT deficiency. Inherited AAT is also commonly associated with liver and biliary disease. Individuals of the heterozygous state have diminished or low normal serum levels of AAT. Approximately 5% to 14% of the adult population is in the heterozygous state and is considered to be at increased risk for the development of emphysema. Homozygous individuals have severe pulmonary and liver disease very early in life. AAT phenotyping is particularly helpful when blood AAT levels are suggestive but not definitive.

Deficiencies of AAT can also be acquired. *Acquired* deficiencies of AAT can occur in patients with protein deficiency syndromes (e.g., malnutrition, liver disease, nephrotic syndrome, and neonatal respiratory distress syndrome).

AAT is also an acute-phase reactant that is elevated in the face of inflammation, infection, or malignancy. It is not specific as to the source of the inflammatory process.

## 38 alpha<sub>1</sub>-antitrypsin

### Interfering factors

- Serum levels of AAT increase during pregnancy or inflammatory diseases.
- ▼ Drugs that may cause *increased* levels include oral contraceptives.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no (verify with laboratory)
- Blood tube commonly used: red
- ▶ If the results show the patient is at risk for developing emphysema, begin patient teaching. Include such factors as avoidance of smoking, infection, and inhaled irritants; proper nutrition; adequate hydration; and education about the disease process of emphysema.

### Abnormal findings

#### ▲ Increased levels

Acute inflammatory disorders  
Chronic inflammatory disorders  
Stress  
Infection  
Thyroid infections

#### ▼ Decreased levels

Early onset of emphysema (in adults)  
Neonatal respiratory distress syndrome  
Cirrhosis (in children)  
Low serum proteins (e.g., nephrotic syndrome, malnutrition, end-stage cancer, protein-losing enteropathy)

notes

## alpha defensin test (Synovasure)

**Type of test** Fluid analysis

### Normal findings

Negative

### Test explanation and related physiology

Infections in prosthetic joints can occur any time after the joint replacement surgery. Infection is a common cause of failed arthroplasties. The diagnosis of periprosthetic joint infection (PJI) remains a serious clinical challenge. Treatment may require prolonged use of parenterally administered antibiotics followed by complex surgery. Therefore it is important that diagnosis is accurate and timely. Diagnosis of PJI is often hampered by administration of antibiotics before or during joint fluid sampling (see p. 118). These antibiotics can inhibit bacterial growth and cause false-negative results.

PJI is usually considered when joint pain worsens after replacement. Joint fluid analysis is not very reliable because of lack of sensitivity and specificity. Advanced joint fluid stains, serology (IgG, IgM), and molecular testing for bacterial DNA sequencing to identify bacteria or the body's response to bacteria have been developed to aid the diagnosis of PJI. Unfortunately, these tests require significant pretest care of the specimen, are time consuming, and too expensive to apply to the population of those potentially affected by PJI. Other testing, such as C-reactive protein (see p. 295) and sed rate (see p. 373), are not specific enough. Imaging (i.e., CT scanning), while more reliable has a low specificity.

Alpha defensin is a protein that is produced by WBCs in joint fluid in response to PJI. An immunoassay to identify this biomarker can accurately identify nearly all incidences of PJI. This testing, however, cannot identify the specific infecting agent or its sensitivity to antibiotics. While synchronously occurring inflammatory disease (e.g., rheumatoid arthritis) can confound other testing for PJI, alpha defensin levels are not affected. A positive alpha defensin test result in the presence of a low C-reactive protein may represent a false-positive result and should be questioned.

### Interfering factors

- Fresh blood in the synovial fluid can alter test results.

## 40 alpha defensin test

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient. Fasting is not required.

- Obtain an informed consent if indicated.

#### During

- After aseptically cleaning and locally anesthetizing the area, a needle is inserted into the affected joint space to remove fluid.
- A few cubic centimeters of synovial fluid is injected into tubes provided by the central laboratory.
- The specimen is sent to the central laboratory.

#### After

**PT** Assess the joint for pain, fever, and swelling. Teach the patient to look for signs of infection at home.

- Apply ice to decrease pain and swelling.

### Abnormal findings

PJI

notes

**alpha-fetoprotein (AFP,  $\alpha_1$ -Fetoprotein)**

**Type of test** Blood

**Normal findings**

Adult: < 40 ng/mL or < 40 mcg/L (SI units)

Child (< 1 year): < 30 ng/mL

(Ranges are stratified by weeks of gestation and vary according to laboratory.)

**Test explanation and related physiology**

Alpha-fetoprotein (AFP) is an oncofetal protein normally produced by the fetal liver and yolk sac. It is the dominant fetal serum protein in the first trimester of life and diminishes to very low levels by the age of 1 year. It is also normally found in very low levels in the adult.

AFP is an effective screening serum marker for fetal body wall defects. The most notable of these are neural tube defects, which can vary from a small myelomeningocele to anencephaly. If a fetus has an open body wall defect, fetal serum AFP leaks out into the amniotic fluid and is picked up by the maternal serum. AFP from fetal sources can normally be detected in the amniotic fluid or the mother's blood after 10 weeks' gestation. Peak levels occur between 16 and 18 weeks' gestation. Maternal serum reflects the changes in amniotic AFP levels. When elevated maternal serum AFP levels are identified, further evaluation with repeat serum AFP levels, amniotic fluid AFP levels, and ultrasound is warranted.

Elevated serum AFP levels in pregnancy may also indicate multiple pregnancy, fetal distress, fetal congenital abnormalities, or intrauterine death. Low AFP levels after correction for age of gestation, maternal weight, race, and presence of diabetes are found in mothers carrying fetuses with trisomy 21 (Down syndrome). See *maternal screen testing* (p. 612) and *nuchal translucency* (p. 685) for other pregnancy screening tests.

AFP is also used as a tumor marker; see p. 194.

**Interfering factors**

- Fetal blood contamination, which may occur during amniocentesis, can cause increased AFP levels.
- Recent administration of radioisotopes can affect values.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no

## 42 alpha-fetoprotein

- Blood tube commonly used: red
- If AFP is to be performed on amniotic fluid, follow the procedure and patient care for amniocentesis (see p. 49).
- Include the gestational age on the laboratory slip.

### Abnormal findings

#### ▲ Increased maternal AFP levels

Neural tube defects  
(e.g., anencephaly,  
encephalocele, spina bifida,  
myelomeningocele)

Abdominal wall defects  
(e.g., gastroschisis or  
omphalocele)

Multiple pregnancy

Threatened abortion

Fetal distress or congenital  
anomalies

Fetal death

#### ▲ Increased nonmaternal AFP levels

Primary hepatocellular cancer (hepatoma)

Germ cell or yolk sac cancer of the ovary

Embryonal cell or germ cell tumor of the testes

Other cancers (e.g., stomach, colon, lung, breast, or  
lymphoma)

Liver cell necrosis (e.g., cirrhosis or hepatitis)

#### ▼ Decreased maternal serum AFP levels

Trisomy 21 (Down  
syndrome)  
Fetal wastage

### notes

**aluminum**

A

**Type of test** Blood**Normal findings**

All ages: 0-6 ng/mL

Dialysis patients of all ages: &lt; 60 ng/mL

**Test explanation and related physiology**

Under normal physiologic conditions, the usual daily dietary intake of aluminum (5-10 mg) is completely excreted by the kidneys. Patients in renal failure (RF) lose the ability to clear aluminum and are at risk for aluminum toxicity. Aluminum-laden dialysis water and aluminum-based phosphate binder gels designed to decrease phosphate accumulation increase the incidence of aluminum toxicity in RF patients. Furthermore, the dialysis process is not highly effective at eliminating aluminum.

If aluminum accumulates, it binds to albumin and is rapidly distributed throughout the body. Aluminum overload leads to accumulation of aluminum in the brain and bone. Brain deposition has been implicated as a cause of dialysis dementia. In bone, aluminum replaces calcium and disrupts normal osteoid formation.

Serum aluminum concentrations are likely to be increased higher than the reference range in patients with metallic joint prosthesis. Serum concentrations greater than 10 ng/mL in a patient with an aluminum-based implant suggest significant prosthesis wear. Chromium and other metals can be detected using similar laboratory techniques.

**Interfering factors**

- Special evacuated blood collection tubes are required for aluminum testing.
- Most of the common evacuated blood collection devices have rubber stoppers that are composed of aluminum silicate. Simple puncture of the rubber stopper for blood collection is sufficient to contaminate the specimen with aluminum.
- Gadolinium- or iodine-containing contrast media that have been administered within 96 hours can alter test results for heavy metals, including aluminum.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no

## 44 aluminum

- Blood tube commonly used: royal blue or tan
- If the blood sample is sent to a central diagnostic laboratory, results will be available in 7 to 10 days.

### Abnormal findings

#### ▲ Increased levels

Aluminum toxicity

notes



**amino acid profiles** (Amino acid screen)

**Type of test** Blood; urine

**Normal findings**

Normal values vary for different amino acids.

**Test explanation and related physiology**

Amino acids are “building blocks” of proteins, hormones, nucleic acids, and pigments. They can act as neurotransmitters, enzymes, and coenzymes. There are eight essential amino acids that must be provided to the body by the diet. The essential amino acids must be transported across the gut and renal tubular lining cells. The metabolism of the essential amino acids is critical to the production of other amino acids, proteins, carbohydrates, and lipids.

When there is a defect in the metabolism or transport of any one of these amino acids, excesses of their precursors or deficiencies of their “end product” amino acid are evident in the blood and/or urine.

Clinical manifestations of these diseases may be precluded if diagnosis is early and if appropriate dietary replacement of missing amino acids is provided. Usually, urine testing for specific amino acids is used to screen for some of these errors in amino acid metabolism and transport. Blood testing is very accurate. Federal law now requires hospitals to test all newborns for inborn errors in metabolism including amino acids. Testing is required for errors in amino acid metabolism such as phenylketonuria (PKU), maple syrup urine disease (MSUD), and homocystinuria (see newborn metabolic screening, p. 642).

A few drops of blood are obtained from the heel of a newborn to fill a few circles on filter paper labeled with names of infant, parent, hospital, and primary physician. The sample is usually obtained on the second or third day of life, after protein-containing feedings (i.e., breast milk or formula) have started.

After a presumptive diagnosis is made, amino acid levels can be detected in blood or amniotic fluid. The genetic defects for many of these diseases are becoming more defined, allowing for even earlier diagnosis to be made in utero.

**Interfering factors**

- The circadian rhythm affects amino acid levels. Levels are usually lowest in the morning and highest by midday.
- Pregnancy is associated with reduced levels of some amino acids.

## 46 amino acid profiles

- ✔ Drugs that may *increase* amino acids include bismuth, heparin, steroids, and sulfonamides.
- ✔ Drugs that may *decrease* some amino acid levels include estrogens and oral contraceptives.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: yes
  - Blood tube commonly used: red
  - Obtain a history of the patient's symptoms.
  - Obtain a pedigree highlighting family members with amino acid disorders.
- PT** A 12-hour fast is generally required before blood collection.

### During

- Usually a 24-hour random urine specimen is required. See inside front cover for Routine Urine Testing.
- Screening is done on a spot urine using the first voided specimen in the morning.

### After

- Generally, genetic counseling is provided before testing and continued after results are obtained.

### Abnormal findings

#### ▲ Increased blood levels

Specific aminoacidopathies  
(e.g., PKU, maple syrup disease)

Specific aminoacidemias  
(e.g., glutaric aciduria)

#### Increased urine levels

Specific aminoacidurias  
(e.g., cystinuria, homocystinuria)

#### ▼ Decreased blood levels

Hartnup disease  
Nephritis  
Nephrotic syndromes

### notes

## ammonia level

**Type of test** Blood

### Normal findings

Adult: 10-80 mcg/dL or 6-47  $\mu\text{mol/L}$  (SI units)

Child: 40-80 mcg/dL

Newborn: 90-150 mcg/dL

### Test explanation and related physiology

Ammonia is used to support the diagnosis of severe liver diseases (fulminant hepatitis or cirrhosis). Ammonia levels are also used in the diagnosis and follow-up of hepatic encephalopathy.

Ammonia is a byproduct of protein catabolism. Most of the ammonia is made by bacteria acting on proteins present in the gut. By way of the portal vein, ammonia goes to the liver where it is normally converted into urea and then secreted by the kidneys. With severe hepatocellular dysfunction, ammonia cannot be catabolized. Furthermore, when portal blood flow to the liver is altered (e.g., in portal hypertension), ammonia cannot reach the liver to be catabolized. Ammonia levels in the blood rise. Plasma ammonia levels do not correlate well with the degree of hepatic encephalopathy. Inherited deficiencies of urea cycle enzymes, inherited metabolic disorders of organic acids, and the dibasic amino acids lysine and ornithine are major causes of high ammonia levels in infants and adults. Finally, impaired renal function diminishes excretion of ammonia, and blood levels rise. High levels of ammonia are often associated with encephalopathy and coma.

### Interfering factors

- Hemolysis increases ammonia levels because the RBCs contain about three times the ammonia content of plasma.
- Muscular exertion can increase ammonia.
- Cigarette smoking can produce significant increases in levels.
- Ammonia levels may be falsely increased if the tourniquet is too tight for a long period.
- ✚ Drugs that may cause *increased* ammonia levels include acetazolamide, alcohol, ammonium chloride, barbiturates, diuretics (e.g., loop, thiazide), narcotics, and parenteral nutrition.
- ✚ Drugs that may cause *decreased* levels include broad-spectrum antibiotics (e.g., neomycin), lactobacillus, lactulose, levodopa, and potassium salts.

## 48 ammonia level

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: green
- Note that some institutions require that the specimen be sent to the laboratory in an iced container.
- Avoid hemolysis and send the specimen promptly to the laboratory.
- Many patients with liver disease have prolonged clotting times.

### Abnormal findings

#### ▲ Increased levels

Primary hepatocellular disease  
Reye syndrome  
Asparagine intoxication  
Portal hypertension  
Severe heart failure with congestive hepatomegaly  
Hemolytic disease of the newborn (erythroblastosis fetalis)  
Gastrointestinal bleeding with mild liver disease  
Gastrointestinal obstruction with mild liver disease  
Hepatic encephalopathy and hepatic coma  
Genetic metabolic disorder of the urea cycle

#### ▼ Decreased levels

Essential or malignant hypertension  
Hyperornithinemia

notes

## amniocentesis (Amniotic fluid analysis)

A

**Type of test** Fluid analysis

### Normal findings

Weeks' gestation	Amniotic fluid volume (mL)
15	450
25	750
30–35	1500
Full term	> 1500

Amniotic fluid appearance: clear; pale to straw yellow

L/S ratio:  $\geq 2:1$

Bilirubin:  $< 0.2$  mg/dL

No chromosomal or genetic abnormalities

Phosphatidylglycerol (PG): positive for PG

Lamellar body count:  $> 30,000$

Alpha-fetoprotein: dependent on gestational age and laboratory technique

Fetal lung maturity (FLM):

Mature:  $< 260$  mPOL

Transitional: 260 to 290 mPOL

Immature:  $> 290$  mPOL

### Test explanation and related physiology

Amniocentesis is performed on pregnant women to gather information about their fetuses. The following can be evaluated by studying the amniotic fluid:

- **Fetal maturity status**, especially pulmonary maturity (when early delivery is preferred). Fetal maturity is determined by analysis of the amniotic fluid in the following manner:
  - a. **Lecithin/sphingomyelin (L/S) ratio.** The L/S ratio is a measure of fetal lung maturity, which is determined by measuring the phospholipids in amniotic fluid. Lecithin is the major constituent of surfactant, an important substance required for alveolar ventilation. If surfactant is insufficient, the alveoli collapse during expiration. This may result in respiratory distress syndrome (RDS), which is a major cause of death in immature babies. In immature fetal lungs, the sphingomyelin concentration in amniotic fluid is higher than the lecithin concentration. At 35 weeks of gestation, whereas the concentration of lecithin rapidly increases, the sphingomyelin concentration decreases. An L/S ratio of 2:1 (3:1 in mothers with diabetes) or greater

is a highly reliable indication that the fetal lungs and therefore the fetus are mature. In such a case, the infant would be unlikely to develop RDS after birth. As the L/S ratio decreases, the risk of RDS increases.

As an alternative to measuring the L/S ratio, the *fetal lung maturity (FLM) test* is based on fluorescence depolarization. This test determines the ratio of surfactant to albumin to evaluate pulmonary maturity.

- b. *Phosphatidylglycerol (PG)*. This is a minor component (~10%) of lung surfactant phospholipids and therefore, alone, is less accurate in measuring pulmonary maturity. However, because PG is almost entirely synthesized by mature lung alveolar cells, it is a good indicator of lung maturity. In healthy pregnant women, PG appears in amniotic fluid after 35 weeks of gestation, and levels gradually increase until term. The simultaneous determination of the L/S ratio and the presence of PG is an excellent method of assessing fetal maturity based on pulmonary surfactant.
- c. *Lamellar body count*. This test to determine fetal maturity is also based on the presence of surfactant. These lamellar bodies represent the storage form of pulmonary surfactant. Lamellar body results are calculated in units of particle density per microliter of amniotic fluid. Some researchers have recommended cutoffs of 30,000/mL and 10,000/mL to predict low and high risks for RDS, respectively. If the count is greater than 30,000, there is a 100% chance that the infant's lungs are mature enough to not experience RDS. If the lamellar body count is less than 10,000, the probability of RDS is high (67%). Values between 10,000 and 30,000/mL represent intermediate risk for RDS.
- d. *Measurement of surfactant activity*. Surfactant activity is a semiquantitative group of tests performed by determining the development and stability of foam when amniotic fluid is shaken in a solution of alcohol. This testing may be called the *tap test*, the *shake test*, or the *foam stability index test*. If a ring of bubbles forms on the surface of the solution, fetal lung maturity is indicated. If no bubbles are present, varying levels of respiratory distress syndrome are indicated.
- e. *Measurement of optical density of amniotic fluid*. The measurement of amniotic fluid at 650 nm can be measured by absorbance. A denser fluid will be associated with greater lung maturity. This testing method is often used as a rapid screening test for fetal lung maturity.

- **Sex of the fetus.** Sons of mothers who are known to be carriers of X-linked recessive traits would have a 50:50 chance of inheritance.
- **Genetic and chromosomal aberrations.** Genetic and chromosomal studies performed on cells aspirated within the amniotic fluid can indicate the gender of the fetus (important in sex-linked diseases such as hemophilia) or the existence of many genetic and chromosomal aberrations (e.g., trisomy 21).
- **Fetal status affected by Rh isoimmunization.** Mothers with Rh isoimmunization may have a series of amniocentesis procedures during the second half of pregnancy to assess the level of bilirubin pigment in the amniotic fluid. The quantity of bilirubin is used to assess the severity of hemolysis in Rh-sensitized pregnancy. Amniocentesis is usually initiated at 24 to 25 weeks' gestation when hemolysis is suspected.
- **Hereditary metabolic disorders,** such as cystic fibrosis.
- **Anatomic abnormalities,** such as neural tube closure defects (myelomeningocele, anencephaly, spina bifida). Increased levels of alpha-fetoprotein (AFP) in the amniotic fluid may indicate a neural crest abnormality (see p. 41). Decreased AFP may be associated with increased risk of trisomy 21.
- **Fetal distress,** detected by meconium staining of the amniotic fluid. This is caused by relaxation of the anal sphincter. In this case the normally colorless or pale, straw-colored amniotic fluid may be tinged with green. Other color changes may also indicate fetal distress. There are, however, more accurate and safer methods of determining fetal stress such as the fetal biophysical profile (see p. 406).
- **Assessment of amniotic fluid for infection.** Amniocentesis is used to obtain fluid for bacterial culture and sensitivity when infection is suspected. This is especially helpful if premature membrane rupture is suspected. Amniotic fluid can also be obtained if viral infections that may affect the fetus are suspected during pregnancy.
- **Assessment for pre-mature rupture of membranes.** Through amniocentesis, a dye can be injected into the amniotic fluid. If this same dye is found in vaginal fluid, rupture of the amniotic membrane is documented. This is sometimes referred to as the *amnio-dye test*. There are, however, more practical tests of vaginal fluid to determine membrane rupture. Most commonly, the pH of the vaginal fluid is determined using a *Nitrazine test* strip. If the test strip turns dark or blue, amniotic fluid is present in the vagina, and membrane rupture is documented. There are now several commercial companies that provide easily performed immunoassays that are able to differentiate normal

## 52 amniocentesis

vaginal discharge from amniotic fluid associated with rupture of membranes. These immunoassays use monoclonal antibodies to identify *placental alpha microglobulin-1* (PAMG-1). The concentration of PAMG-1 in vaginal fluid is thousands of times higher in amniotic fluid than it is in normal vaginal secretions, thus allowing differentiation of the two.

The timing of the amniocentesis varies according to the clinical circumstances. With advanced maternal age and if chromosomal or genetic aberrations are suspected, the test should be done early enough (at 14 to 16 weeks of gestation; at least 150 mL of fluid exists at this time) to allow a safe abortion. If information on fetal maturity is sought, performing the study during or after the 35th week of gestation is best.

### Contraindications

- Patients with abruptio placentae
- Patients with placenta previa
- Patients with a history of premature labor (before 34 weeks of gestation unless the patient is receiving antilabor medication)
- Patients with an incompetent cervix or cervical insufficiency
- Patients with anhydramnios
- Patients with suspected premature labor

### Potential complications

- Miscarriage
- Fetal injury
- Leak of amniotic fluid
- Infection (amnionitis)
- Abortion
- Premature labor
- Maternal hemorrhage
- Maternal Rh isoimmunization
- Amniotic fluid embolism
- Abruptio placentae
- Inadvertent damage to the bladder or intestines

### Interfering factors

- Fetal blood contamination can cause false AFP elevations.
- Hemolysis of the specimen can alter results.
- Contamination of the specimen with meconium or blood may give inaccurate L/S ratios.

### Procedure and patient care

#### Before


- PT** Explain the procedure to the patient. Allay any fears and allow the patient to verbalize her concerns.



- Obtain an informed consent.
- PT** Tell the patient that no food or fluid is restricted.
- Evaluate the mother's blood pressure and the fetal heart rate.
  - Follow instructions regarding emptying the bladder, which depend on gestational age. Before 20 weeks of gestation, the bladder may be kept full to support the uterus. After 20 weeks, the bladder may be emptied to minimize the chance of puncture.
  - Note that the placenta is localized before the study by ultrasound to permit selection of a site that will avoid placental puncture.

### During

- Place the patient in the supine position.
- Note the following procedural steps:
  1. The skin overlying the chosen site is prepared and usually anesthetized locally.
  2. A needle with a stylet is inserted through the midabdominal wall and is directed at an angle toward the middle of the uterine cavity (Figure 2).

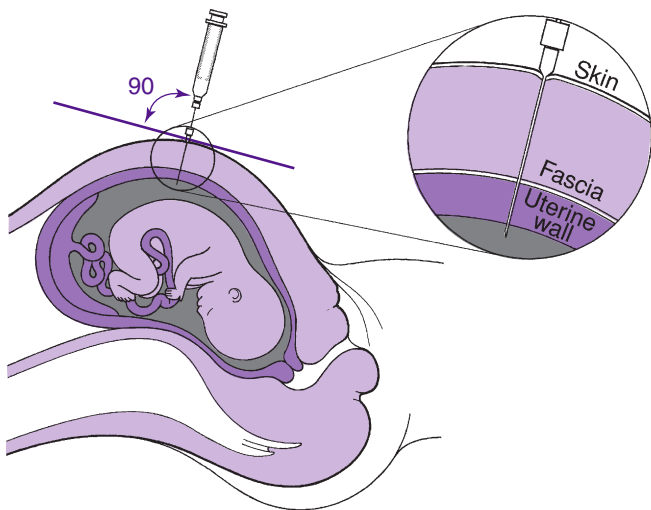

**FIG. 2** Amniocentesis. Ultrasound scanning is usually used to determine the placental site and to locate a pocket of amniotic fluid. The needle is then inserted. Three levels of resistance are felt as the needle penetrates the skin, fascia, and uterine wall. When the needle is placed within the uterine cavity, amniotic fluid is withdrawn.

## 54 amniocentesis

3. The stylet is then removed and a sterile plastic syringe attached.
  4. After 5 to 10 mL of amniotic fluid is withdrawn, the needle is removed.
  5. The specimen is placed in a light-resistant container to prevent breakdown of bilirubin.
  6. The site is covered with an adhesive bandage.
  7. If the amniotic fluid is bloody, the physician must determine whether the blood is maternal or fetal in origin. The Kleihauer–Betke stain will stain fetal cells pink. Meconium in the fluid is usually associated with a compromised fetus.
  8. Amniotic fluid volume is calculated.
- Note that this procedure is performed by a physician and takes approximately 20 to 30 minutes.
- PT** Tell the patient that the discomfort associated with amniocentesis is usually described as a mild uterine cramping that occurs when the needle contacts the uterus. Some women may complain of a pulling sensation as the amniotic fluid is withdrawn.
- Remember that many women are extremely anxious during this procedure.

### After

- Place amniotic fluid in a sterile, siliconized glass container and transport it to a special chemistry laboratory for analysis. Sometimes the specimen may be sent by airmail to another commercial laboratory.
- PT** Inform the patient that the results of this study are usually not available for more than 1 week.
- For women who have Rh-negative blood, administer RhoGAM because of the risk of isoimmunization from the fetal blood.
  - Assess the fetal heart rate after the test to detect any ill effects related to the procedure. Compare this value with the preprocedure baseline value.
- PT** If the patient felt dizzy or nauseated during the procedure, instruct her to lie on her left side for several minutes before leaving the examining room.
- Observe the puncture site for bleeding or other drainage.
- PT** Instruct the patient to call her physician if she has any amniotic fluid loss, bleeding, temperature elevation, abdominal pain, abdominal cramping, fetal hyperactivity, or unusual fetal lethargy.

## Abnormal findings

Hemolytic disease of the newborn

Rh isoimmunization

Neural tube closure defects (e.g., myelomeningocele, anencephaly, spina bifida)

Abdominal wall closure defects (e.g., gastroschisis, omphalocele)

Sacrococcygeal teratoma

Meconium staining

Immature fetal lungs

Hereditary metabolic disorders (e.g., cystic fibrosis, Tay-Sachs disease, galactosemia)

Genetic or chromosomal aberrations (e.g., sickle cell anemia, thalassemia, trisomy 21 [Down syndrome])

Sex-linked disorders (e.g., hemophilia)

Polyhydramnios

Oligohydramnios

## notes

## amylase

**Type of test** Blood; urine

### Normal findings

#### Blood

Adult: 60-120 Somogyi units/dL or 30-220 units/L (SI units)

Values may be slightly increased during normal pregnancy and in the elderly.

Newborn: 6-65 units/L

#### Urine (24-hour)

Up to 5000 Somogyi units/24 hr or 6.5-48.1 units/hr (SI units)

### Possible critical values

**Blood:** More than three times the upper limit of normal (depending on the method)

### Test explanation and related physiology

*Serum amylase* is an easily and rapidly performed test that is commonly used to diagnose and monitor the treatment of pancreatitis or obstruction of the pancreatic duct flow (as a result of pancreatic carcinoma). Blood vessels draining the free peritoneum and absorbing the lymph pick up the excess amylase. An abnormal rise in the serum level of amylase occurs within 12 hours of the onset of disease. Because amylase is rapidly cleared by the kidneys, serum levels return to normal 48 to 72 hours after the initial insult. Persistent pancreatitis, duct obstruction, or pancreatic duct leak will cause persistent elevated amylase levels.

Although serum amylase is a sensitive test for pancreatic disorders, it is not specific. Other nonpancreatic diseases can cause elevated amylase levels in the serum. For example, in a bowel perforation, intraluminal amylase leaks into the free peritoneum and is picked up by the peritoneal blood vessels. Also, a penetrating peptic ulcer into the pancreas will cause elevated amylase levels. Duodenal obstruction can be associated with less significant elevations in amylase. Because salivary glands contain amylase, elevations can be expected in patients with parotiditis (mumps). Amylase isoenzyme testing can differentiate pancreatic from salivary hyperamylasemia.

*Urine amylase* levels rise after the blood levels. Several days after the onset of the disease process, serum amylase levels may be normal, but urine amylase levels are significantly elevated. Urine amylase is particularly useful in detecting pancreatitis late in the disease course.

As with serum amylase, urine amylase is sensitive but not specific for pancreatic disorders. A comparison of the renal clearance ratio of amylase with creatinine provides more specific diagnostic information than either the urine amylase level or the serum amylase level alone. When the *amylase/creatinine clearance* ratio is 5% or more, the diagnosis of pancreatitis can be made with certainty.

### Interfering factors

- Serum lipemia may falsely decrease amylase levels.
- ✚ IV dextrose solutions can cause a false-negative result.
- ✚ Drugs that may cause *increased* serum amylase levels include aminosalicic acid, aspirin, azathioprine, corticosteroids, dexamethasone, ethyl alcohol, glucocorticoids, iodine-containing contrast media, loop diuretics, methyldopa, narcotic analgesics, oral contraceptives, and prednisone.
- ✚ Drugs that may cause *decreased* levels include citrates, glucose, and oxalates.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- See inside front cover for Routine Urine Testing.

### Abnormal findings

#### ▲ Increased levels

Acute or chronic relapsing pancreatitis  
 Penetrating or perforated peptic ulcer  
 Necrotic or perforated bowel  
 Acute cholecystitis  
 Parotiditis (mumps)  
 Ectopic pregnancy  
 Pulmonary infarction  
 Diabetic ketoacidosis  
 Duodenal obstruction  
 Osteogenic sarcoma  
 Cryoglobulinemia  
 Rheumatoid diseases

**amyloid beta protein precursor, soluble** (sBPP, APP)

**Type of test** Cerebrospinal fluid (CSF) analysis

**Normal findings**

> 450 units/L

**Test explanation and related physiology**

This test is used to help diagnose Alzheimer disease (AD) and other forms of senile dementia. Amyloid beta protein is formed by cleavage of amyloid precursor protein (APP). There are many isoforms of sBPP depending on how many amino acids are cleaved into the form (30-51). The 42-amino-acid peptide ( $AB_{42}$ ) and  $AB_{40}$  are most significant to the diagnosis of AD. These beta amyloid proteins have been shown to be neurotrophic and neuroprotective. Beta amyloid is deposited on the brain in the form of plaques in patients with AD. As a result of this deposition, levels of beta amyloid are *decreased* in the cerebrospinal fluid of patients with AD and other forms of dementia. The  $AB_{42} / AB_{40}$  ratio in CSF may be more accurate in the diagnosis of AD. This test is often performed as a diagnostic panel that includes Apolipoprotein E (see p. 101) testing.

Ongoing research has also focused on using CSF levels of *tau protein* as another biochemical marker for AD.

Recently, PET scanning with amyloid imaging has shown promise for the diagnosis of AD. Pittsburgh Agent B (PIB) appears to reliably detect brain amyloid caused by the accumulation of A beta 42 within plaques. Studies so far have revealed high levels of amyloid retention in the brain at prodromal stages of AD and the possibility of discriminating AD from other dementia disorders by scanning with PIB. Because amyloid accumulation is one of the earliest signs of AD, early diagnosis may be facilitated by identifying amyloid early in the disease progression, perhaps before symptoms emerge.

**Procedure and patient care****Before**

**pt** Explain the procedure to the patient.

- Refer to the instructions for a lumbar puncture and CSF examination (see p. 576).

**During**

- Collect a CSF specimen as per the lumbar puncture discussion.

### After

- Follow the postprocedure guidelines after a lumbar puncture.

### Abnormal findings

- ▼ **Decreased levels**  
Alzheimer disease

notes

## angiotensin

**Type of test** Blood

### Normal findings

Angiotensin I:  $\leq 25$  pg/mL

Angiotensin II: 10-60 pg/mL

### Test explanation and related physiology

Renin (see p. 783) is an enzyme that is released by the juxtaglomerular apparatus of the kidneys. Its release is stimulated by hypokalemia, hyponatremia, decreased renal blood perfusion, or hypovolemia. Renin stimulates the release of angiotensinogen. Angiotensin-converting enzyme (ACE) (see p. 61) metabolizes angiotensinogen to angiotensin I and subsequently to angiotensin II and III. Angiotensin then stimulates the release of catecholamines, antidiuretic hormone, ACTH, oxytocin, and aldosterone. Angiotensin is also a vasoconstrictor. Angiotensin is used to identify renovascular sources of hypertension. It can be measured as angiotensin I or as angiotensin II.

### Interfering factors

See renin, p. 783.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender
- PT** Instruct the patient to maintain a normal diet with a restricted amount of sodium ( $\sim 3$  g/day) for 3 days before the test.
- PT** Instruct the patient to check with a healthcare provider about discontinuing any medications that may interrupt renin activity.
- Record the patient's position, dietary status, and time of day on the laboratory slip.
- Place the tube of blood on ice and immediately send it to the laboratory.

### Abnormal findings

#### ▲ Increased levels

Essential hypertension  
Malignant hypertension  
Renovascular  
hypertension

#### ▼ Decreased levels

Primary hyperaldosteronism  
Steroid therapy  
Congenital adrenal  
hyperplasia



## angiotensin-converting enzyme (ACE, Serum angiotensin-converting enzyme [SACE])

**Type of test** Blood

### Normal findings

8-53 U/L

### Test explanation and related physiology

ACE is used to detect and monitor the clinical course of sarcoidosis. It is also used to differentiate sarcoidosis from other granulomatous diseases and to differentiate between active and dormant sarcoid disease.

Elevated ACE levels are found in a high percentage of patients with sarcoidosis. This test is used primarily in patients with sarcoidosis to evaluate the severity of disease and the response to therapy. Levels are especially high with active pulmonary sarcoidosis and can be normal with inactive (dormant) sarcoidosis. Elevated ACE levels also occur in conditions other than sarcoidosis, including Gaucher disease (a rare familial lysosomal disorder of fat metabolism), leprosy, alcoholic cirrhosis, active histoplasmosis, tuberculosis, Hodgkin disease, myeloma, scleroderma, pulmonary embolism, and idiopathic pulmonary fibrosis. ACE is elevated in the CSF of patients with neurosarcoidosis.

### Interfering factors

- Patients younger than 20 years of age normally have very high ACE levels.
- Hemolysis or hyperlipidemia may falsely decrease ACE levels.
- ✚ Drugs that may cause *decreased* ACE levels include ACE inhibitor antihypertensives and steroids.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Note on the laboratory slip if the patient is taking steroids.

### Abnormal findings

#### ▲ Increased levels

Sarcoidosis  
 Gaucher disease  
 Tuberculosis  
 Leprosy  
 Alcoholic cirrhosis

## 62 angiotensin-converting enzyme

### ▲ Increased levels

- Active histoplasmosis
- Hodgkin disease
- Myeloma
- Idiopathic pulmonary fibrosis
- Diabetes mellitus
- Primary biliary cirrhosis
- Amyloidosis
- Hyperthyroidism
- Scleroderma
- Pulmonary embolism

notes

**anion gap (AG, R factor)**

**Type of test** Blood

**Normal findings**

$16 \pm 4$  mEq/L (if potassium is used in the calculation)

$12 \pm 4$  mEq/L (if potassium is not used in the calculation)

**Test explanation and related physiology**

The anion gap (AG) is the difference between the cations and the anions in the extracellular space that is routinely calculated in the laboratory (i.e.,  $AG = [\text{sodium} + \text{potassium}] - [\text{chloride} + \text{bicarbonate}]$ ). In some laboratories, the potassium is not measured because the level of potassium in acid-base abnormalities varies. The normal value of the AG is adjusted downward if potassium is eliminated from the equation. The AG, although not real physiologically, is created by the small amounts of anions in the blood (e.g., lactate, phosphates, sulfates, organic anions, and proteins) that are not measured.

This calculation is most often helpful in identifying the cause of metabolic acidosis. As such acids as lactic acid or ketoacids accumulate in the bloodstream, bicarbonate neutralizes them to maintain a normal pH within the blood. Mathematically, when bicarbonate decreases, the AG increases. In general, most metabolic acidotic states (excluding some types of renal tubular acidosis) are associated with an increased AG. The higher the gap is above normal, the more likely that the metabolic acidotic state is associated with the AG. Proteins can have a significant effect on AG. As albumin (usually negatively charged) increases, AG will increase.

A decreased AG is very rare but can occur when there is an increase in unmeasured (calcium or magnesium) cations. A reduction in anionic proteins (nephrotic syndrome) will also decrease AG. For example, a 1 g/dL drop in serum protein is associated with a 2.5 mEq/L drop in AG. Because the anion proteins are lost, the  $\text{HCO}_3^-$  increases to maintain electrical neutrality. An increase in cationic proteins (some immunoglobulins) will also decrease AG. Except for hypoproteinemia, conditions that cause a reduced or negative AG are relatively rare compared with those associated with an elevated AG.

**Interfering factors**

- Hyperlipidemia may cause undermeasurement of sodium and falsely decrease AG.

## 64 anion gap

- Normal values of AG vary according to different normal values for electrolytes, depending on laboratory methods of measurement.
- ▼ Drugs that *increase* AG are many. Examples include carbonic anhydrase inhibitors (e.g., acetazolamide), ethanol, methanol, and salicylate.
- ▼ Drugs that *decrease* AG are also many. Examples include acetazolamide, lithium, spironolactone, and sulindac.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or green
- If the patient is receiving an IV infusion, obtain the blood from the opposite arm.
- The sodium, potassium, chloride, and bicarbonate levels are determined by an automated multichannel analyzer. The AG is then calculated as indicated in the test explanation section.

### Abnormal findings

#### ▲ Increased levels

Lactic acidosis  
Diabetic ketoacidosis  
Alcoholic ketoacidosis  
Starvation  
Renal failure  
Renal tubular acidosis  
Increased gastrointestinal losses of bicarbonate (e.g., diarrhea or fistulae)  
Hypoaldosteronism

#### ▼ Decreased levels

Excess alkali ingestion  
Multiple myeloma  
Chronic vomiting or gastric suction  
Hyperaldosteronism  
Hypoproteinemia  
Lithium toxicity  
Bromide (cough syrup) toxicity

### notes

## anticardiolipin antibodies (aCL antibodies, ACA, Antiphospholipid antibodies, Lupus anticoagulant)

**Type of test** Blood

### Normal findings

Negative:

- < 23 GPL (IgG phospholipid units)
- < 11 MPL (IgM phospholipid units)

### Test explanation and related physiology

Anticardiolipin antibodies (immunoglobulins G and M to cardiolipin) are antiphospholipid autoantibodies that attach to phospholipids on cell membranes and can interfere with the coagulation system. *Antiphospholipid autoantibodies* include *anticardiolipin antibodies* and the *lupus anticoagulant antibody*. Phospholipid antibodies occur in patients with a variety of clinical signs and symptoms, notably thrombosis (arterial or venous), pregnancy morbidity (unexplained fetal death, premature birth, severe preeclampsia, or placental insufficiency), unexplained cutaneous circulation disturbances (levido reticularis or pyoderma gangrenosum), thrombocytopenia or hemolytic anemia, and nonbacterial thrombotic endocarditis. Phospholipid antibodies and lupus anticoagulants are found with increased frequency in patients with systemic rheumatic diseases, especially lupus erythematosus. The term *antiphospholipid syndrome (APS)* or *Hughes syndrome* is used to describe the triad of thrombosis, recurrent fetal loss, and thrombocytopenia accompanied by phospholipid antibodies or a lupus anticoagulant. These antibodies may be considered normal in an elderly person.

### Interfering factors

- Patients who have or had syphilis infections can have a false-positive result.
- Transient antibodies can occur in patients with infections, AIDS, inflammation, autoimmune diseases, or cancer.
- ✚ False-positive results have been seen in patients who take such medications as chlorpromazine, hydralazine, penicillin, phenytoin, procainamide, and quinidine.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: Verify with laboratory.

## 66 anticardiolipin antibodies

### Abnormal findings

#### ▲ Increased levels

- Systemic lupus erythematosus
- Thrombosis
- Thrombocytopenia
- Recurrent fetal loss
- Syphilis
- Acute infection
- Elderly persons

#### notes

**antacentromere antibody test** (Centromere antibody)

**Type of test** Blood

**Normal findings**

Negative (if positive, serum will be titrated)

Weak positive: positive at the screening titer (1:40 for human epithelial type 2 cells [HEp-2 cells]) (1:20 for kidney cells)

Moderately positive: one dilution higher than screening titer

Strong positive: two dilutions higher than screening titer

**Test explanation and related physiology**

A centromere is the region of the chromosome referred to as the *primary constriction* that divides the chromosome into *arms*. During cell division, the centromere exists in the *pole* of the mitotic spindle.

Anticentromere antibodies are a form of *antinuclear antibodies*. They are found in a very high percentage of patients with CREST syndrome, a variant of scleroderma. CREST syndrome is characterized by calcinosis, Raynaud phenomenon, esophageal dysfunction, sclerodactyly, and telangiectasia. Anticentromere antibodies, on the contrary, are present in only a small minority of patients with scleroderma, a disease that is difficult to differentiate from CREST syndrome. No correlation exists between antibody titer and severity of CREST syndrome.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**Abnormal findings****Positive results**

CREST syndrome

notes

## antichromatin antibody test (Antinucleosome antibody [anti-NCS], Antihistone antibody test [anti-HST, AHA])

**Type of test** Blood

### Normal findings

#### Antinucleosome antibodies

No antibodies present in < 1:20 dilution

#### Antihistone antibody

None detected: < 1.0 units

Inconclusive: 1.0-1.5 units

Positive: 1.6-2.5 units

Strong positive: > 2.5 units

### Test explanation and related physiology

There are several chromatin antinuclear antibodies associated with autoimmune diseases. Nucleosome (NCS) represents the main autoantigen-immunogen in systemic lupus erythematosus (SLE), and these specific antibodies are an important marker of the disease activity. *Antinucleosome (anti-NCS; antichromatin)* antibodies play a key role in the pathogenesis of SLE. Nearly all patients with SLE have anti-NCS antibodies. Anti-NCS is also an antinuclear antibody (see p. 86). Anti-NCS has a sensitivity of 100% and specificity of 97% for SLE diagnosis. Anti-NCS antibodies show the highest correlation with disease activity. Anti-NCS antibodies also show strong association with renal damage (glomerulonephritis and proteinuria) associated with SLE. Anti-NCS autoantibodies are more prevalent than anti-DNA in SLE patients.

Histone antibodies are present in 20% to 55% of idiopathic SLE and 80% to 95% of drug-induced lupus erythematosus. They occur in fewer than 20% of other types of connective tissue diseases. This antibody is particularly helpful in identifying patients with drug-induced lupus erythematosus from drugs such as procainamide, quinidine, penicillamine, hydralazine, methyldopa, isoniazid, and acebutolol. There are several subtypes of *antihistone antibodies (AHAs)*. In drug-induced lupus erythematosus, a specific AHA (anti-[(H2A-H2B)-DNA] IgG) is produced, whereas in most of the other associated diseases (rheumatoid arthritis, juvenile rheumatoid arthritis, primary biliary cirrhosis, autoimmune hepatitis, and dermatomyositis/polymyositis), the AHAs are of other varying specificities.



**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or gold

**Abnormal findings****▲ Increased levels**

- Systemic lupus erythematosus
- Drug-induced lupus erythematosus
- Other autoimmune diseases

**notes**

## anticyclic citrullinated peptide antibody (Cyclic citrullinated peptide antibody, CCP IgG anti-CCP) and anti-mutated citrullinated vimentin (anti-MCV)

**Type of test** Blood

### Normal findings

< 20 units/mL

### Test explanation and related physiology

Anti-CCP is known more formally as *anticyclic citrullinated peptide antibody*. The citrullinated peptide antigen is formed by the intermediary conversion of the amino acid ornithine to arginine. Anti-CCP appears early in the course of rheumatoid arthritis (RA). When the citrulline antibody is detected in a patient's blood, there is a high likelihood that the patient has RA. Anti-CCP is therefore useful in the diagnosis of patients with unexplained joint inflammation, especially when the traditional blood test, rheumatoid factor (RF; see p. 790), is negative or lower than 50 units/mL.

Many patients with early RA may not have elevation of RF, making the diagnosis difficult in the initial stage. The diagnosis of RA can be made, even if RF is negative, if the anti-CCP is elevated. This is particularly important because aggressive treatment in the early stages of RA prevents progression of joint damage. Anti-CCP may occur years before any clinical onset of arthritis or significant elevation of anti-CCP. At a cutoff of 5 units/mL, the sensitivity and specificity of anti-CCP for RA are 67.5% and 99.3%, respectively. RF has a sensitivity of 66.3% and a lower specificity (82.1%) than anti-CCP. When the two antibodies are used together, the specificity for diagnosing RA is 99.1%.

The presence of anti-CCP in RA indicates a more aggressive and destructive form of the disease. It is also a marker for disease progression.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Rheumatoid arthritis

### notes

**antidiuretic hormone (ADH, Vasopressin)**

**Type of test** Blood

**Normal findings**

ADH: 1-5 pg/mL or 1-5 ng/L (SI units)

ADH suppression test (water load test):

65% of water load excreted in 4 hours

80% of water load excreted in 5 hours

Urine osmolality (in second hour)  $\leq$  100 mmol/kg

Urine to serum (U/S) osmolality ratio  $>$  100

Urine specific gravity  $<$  1.003

**Test explanation and related physiology**

ADH, also known as *vasopressin*, is formed by the hypothalamus and is stored in the posterior pituitary gland. It controls the amount of water resorbed by the kidneys. ADH release is stimulated by an increase in serum osmolality or a decrease in intravascular blood volume. Physical stress, surgery, and even high levels of anxiety may also stimulate ADH release. With a release of ADH, more water is resorbed from the kidneys. This increases the amount of free water in the bloodstream and causes a very concentrated urine. With low ADH levels, water is allowed to be excreted, thereby producing hemoconcentration and a more dilute urine.

Diabetes insipidus (DI) results when ADH secretion is inadequate or when the kidney is unresponsive to ADH stimulation. Inadequate ADH secretion is usually associated with central neurologic abnormalities (neurogenic DI), such as trauma, tumor, inflammation of the brain (hypothalamus), or surgical ablation of the pituitary gland. Patients with DI excrete large volumes of free water within a dilute urine. Their blood is hemoconcentrated, causing them to have a strong thirst.

Primary renal diseases may make the renal collecting system less sensitive to ADH stimulation (nephrogenic DI). Again, in this instance, a dilute urine created by excretion of high volumes of free water may occur. To differentiate neurogenic DI from nephrogenic DI, a *water deprivation test (ADH stimulation test)* is performed. During this test, water intake is restricted, and urine osmolality is measured before and after vasopressin is administered. In neurogenic DI, there is no rise in urine osmolality with water restriction, but there is a rise after vasopressin administration. In nephrogenic DI, there is no rise in urine osmolality after water deprivation or vasopressin administration. The diagnosis

## 72 antidiuretic hormone

indicated by this test can be corroborated by a serum ADH level. In neurogenic DI, ADH levels are low. In nephrogenic DI, ADH levels are high.

High serum ADH levels are also associated with the syndrome of inappropriate antidiuretic hormone (SIADH) secretion. In response to the inappropriately high level of ADH secretion, water is resorbed by the kidneys greatly in excess of normal amounts. Thus the patient becomes very hemodiluted and the urine concentrated. Blood levels of important serum ions diminish, causing severe neurologic, cardiac, and metabolic alterations. SIADH can be associated with pulmonary diseases (e.g., tuberculosis, bacterial pneumonia), severe stress (e.g., surgery, trauma), CNS tumor, or infection. Ectopic secretion of ADH from neoplasm (paraneoplastic syndrome) can cause SIADH. The most common tumors associated with SIADH include carcinomas of the lung and thymus, lymphomas, leukemia, and carcinomas of the pancreas, urologic tract, and intestine. Patients with myxedema or Addison disease also can experience SIADH. Some drugs are also known to cause SIADH (see Interfering factors).

The *water load test (ADH suppression test)* is used to differentiate SIADH from other causes of hyponatremia or edematous states. Usually this test is done concomitant with measurements of urine and serum osmolality. Patients with SIADH will excrete none or very little of the water load. Furthermore, their urine osmolality will never be less than 100, and the urine/serum ratio is greater than 100. Patients with other hyponatremia, edematous states, or chronic renal diseases will excrete up to 80% of the water load and will develop midrange osmolality results.

### Interfering factors

- Patients with dehydration, hypovolemia, and stress may have increased ADH levels.
- Patients with overhydration, decreased serum osmolality, and hypervolemia may have decreased ADH levels.
- Use of a glass syringe or collection tube causes degradation of ADH.
- ✎ Drugs that *increase* ADH levels and may cause SIADH include acetaminophen, barbiturates, carbamazepine, cholinergic agents, cyclophosphamide, some diuretics (e.g., thiazides), estrogen, narcotics, nicotine, oral hypoglycemic agents (particularly sulfonylureas), and tricyclic or SSRI antidepressants.
- ✎ Drugs that *decrease* ADH levels include alcohol, beta-adrenergic agents, morphine antagonists, and phenytoin.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red
- Evaluate the patient for high levels of physical or emotional stress.
- Collect a venous blood sample in a *plastic* red-top tube with the patient in the sitting or recumbent position.
- The *water load* test necessitates a baseline serum sodium before water administration. Urine is later collected for specific gravity and osmolality. Blood is collected for osmolality.

## Abnormal findings

### ▲ Increased levels

Syndrome of inappropriate antidiuretic hormone (SIADH)

Nephrogenic diabetes insipidus caused by primary renal diseases

Postoperative days 1-3

Severe physical stress (e.g., trauma, pain, prolonged mechanical ventilation)

Hypovolemia

Dehydration

Acute porphyria

### ▼ Decreased levels

Neurogenic (or central) diabetes insipidus

Surgical ablation of the pituitary gland

Hypervolemia

Decreased serum osmolality

## notes

**anti-DNA antibody test** (Anti-deoxyribonucleic acid antibodies, Antibody to double-stranded DNA, Anti-double-stranded DNA, Anti-ds-DNA, DNA antibody, Native double-stranded DNA)

**Type of test** Blood

### Normal findings

Negative: < 5 international units/mL

Intermediate: 5-9 international units/mL

Positive:  $\geq$  10 international units/mL

### Test explanation and related physiology

The anti-DNA test is useful for the diagnosis and follow-up of systemic lupus erythematosus (SLE). This antibody is found in approximately 65% to 80% of patients with active SLE and rarely in other diseases. The anti-DNA titer decreases with successful therapy and increases with an exacerbation of SLE, especially with the onset of lupus glomerulonephritis.

The anti-DNA IgG antibody is a subtype of the *antinuclear antibodies (ANAs)* (see p. 86). If the ANAs are negative, there is no reason to test for anti-DNA antibodies. There are two types of anti-DNA antibodies. The first and most popular is the antibody against double-stranded DNA (anti-ds-DNA). The second type is the antibody against single-stranded DNA (anti-ss-DNA), which is less sensitive and specific for SLE but is positive in other autoimmune diseases.

### Interfering factors

- A radioactive scan performed within 1 week before the test may alter the test results.
- ✦ Drugs that may cause *increased* levels include hydralazine and procainamide.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Collagen vascular disease (e.g., systemic lupus erythematosus)  
 Chronic hepatitis  
 Infectious mononucleosis  
 Biliary cirrhosis

**antiextractable nuclear antigens** (Anti-ENAs, Antibodies to extractable nuclear antigens, Antihistidyl transfer synthase [anti-Jo-1], Antiribonucleoprotein [anti-RNP], Anti-Smith [anti-SM])

**Type of test** Blood

### Normal findings

Negative

### Test explanation and related physiology

The anti-ENAs are used to assist in the diagnosis of systemic lupus erythematosus (SLE) and mixed connective tissue disease (MCTD) and to eliminate other rheumatoid diseases.

Anti-ENAs are a type of *antinuclear antibody* to certain nuclear antigens that consist of RNA and protein. The ENA antigen is sometimes referred to as *saline-extracted antigen*. The most common ENAs are *Smith (SM)* and *ribonucleoprotein (RNP)*.

The *antinuclear Smith (anti-SM)* antibody is present in about 30% of patients with SLE and in about 8% of patients with MCTD. However, it is not present in patients with most other rheumatoid-collagen diseases.

The *antinuclear ribonucleoprotein (anti-RNP)* antibody is reported in nearly 100% of patients with MCTD and in about 25% of patients with SLE, discoid lupus, and progressive systemic sclerosis (scleroderma). In high titer, anti-RNP is suggestive of MCTD.

The *antihistidyl transfer synthase (anti-Jo-1)* antibodies occur in patients with autoimmune interstitial pulmonary fibrosis and in a minority of patients with aggressive autoimmune myositis.

There are two other antibodies to ENAs. Anti-SS-A and anti-SS-B are described on p. 95 and are used mainly in the diagnostic evaluation of Sjögren syndrome.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

- ▲ **Increased anti-SM antibodies**  
Systemic lupus erythematosus

## 76 antiextractable nuclear antigens

### ▲ Increased anti-RNP antibodies

Mixed connective tissue disease  
Systemic lupus erythematosus  
Discoid lupus scleroderma

### ▲ Increased anti-Jo-1 antibodies

Pulmonary fibrosis  
Autoimmune myositis

notes



**antifactor Xa (Anti-Xa)**

**Type of test** Blood

**Normal findings****Adults and children > 8 weeks**

Therapeutic ranges of heparin:

LMWH: 0.5-1.2 IU/mL

UFH: 0.3-0.7 IU/mL

Prophylactic ranges of heparin:

LMWH: 0.2-0.5 IU/mL

UFH: 0.1-0.4 IU/mL

**Children <8 weeks**

Therapeutic ranges:

Standard heparin (UFH): 0.3 to 0.7 IU/mL

LMWH: 0.5 to 1 IU/mL (specimen drawn 4-6 hours after subcutaneous injection)

Prophylactic range:

LMWH: 0.1 to 0.3 IU/mL

**Test explanation and related physiology**

Through its action on the protein antithrombin, heparin interferes with the clotting process by accelerating the inhibition of coagulation factors Xa and IIa. The plasma Anti-Xa assay (see secondary hemostasis, [Figure 10](#), p. 251) is used to monitor patients on low molecular weight heparins (LMWHs) (e.g., enoxaparin or rivaroxaban) or unfractionated standard heparin (UFH). Although UFH is commonly monitored by means of the APTT (see p. 681), the APTT can underestimate the degree of heparin anticoagulation in patients with a high factor VIII level or suspected heparin resistance (also evaluated by APC; see p. 397). In these situations, the measurement of a plasma anti-Xa level may provide a more accurate assessment of anticoagulation. Reference ranges for anti-Xa levels vary depending on the laboratory performing the test. When a person is not taking heparin, anti-Xa concentrations should be zero or undetectable.

When used to monitor LMWHs, anti-Xa levels are usually ordered as a “peak” test. They are collected 3 to 4 hours after a LMWH dose is given when the blood level is expected to be the highest. Random and “trough” anti-Xa tests may also be ordered when there are concerns that a LMWH may be accumulating in patients with renal failure.

## 78 antifactor Xa

As with all heparins, protamine sulfate can inhibit the drug action. The following formula is useful for calculating the dose of protamine sulfate required to neutralize UFH based upon the anti-Xa level:

$$\text{Anti Xa (IU/mL)} \times \text{Plasma volume (mL/kg)} = \text{dose of Protamine Sulfate (mg)}$$

### Interfering factors

- ✓ Accurate testing requires indication of which LMWH is being administered. Testing methods vary based on the LMWH used.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: blue
- If the patient is receiving heparin by intermittent injection, plan to draw the blood specimen 30 minutes to 1 hour before the next dose of heparin.
- If the patient is receiving continuous heparin, draw the blood at any time.
- When it is used to monitor LMWHs, the specimen is collected 3 to 4 hours after a LMWH dose is administered.
- Apply pressure to the venipuncture site. Remember, if the patient is receiving anticoagulants or has coagulopathies, the bleeding time will be increased.
- Assess the patient to detect possible bleeding. Check for blood in the urine and all other excretions and assess the patient for bruises, petechiae, and low back pain.
- If severe bleeding occurs, note that the anticoagulant effect of heparin can be reversed by parenteral administration of protamine sulfate.

### Abnormal findings

#### ▲ Increased levels

Heparin administration  
Heparin resistance  
Renal failure  
Hemophilia  
Lupus anticoagulant

notes

## antiglomerular basement membrane antibodies

(Anti-GBM antibody, AGBM, Glomerular basement antibody, Goodpasture antibody)

**Type of test** Blood; microscopic examination of tissue

### Normal findings

#### Tissue

Negative: no immunofluorescence noted on the renal or lung tissue basement membrane

#### Blood: EIA (enzyme immunoassay)

Negative: < 20 units

Borderline: 20 to 100 units

Positive: > 100 units

### Test explanation and related physiology

This test is used to detect the presence of circulating glomerular basement membrane (GBM) antibodies commonly present in autoimmune-induced nephritis (Goodpasture syndrome).

Goodpasture syndrome is an autoimmune disease characterized by the presence of antibodies circulating against antigens in the basement membrane of the renal glomerular and the pulmonary alveoli. These immune complexes activate the complement system and thereby cause tissue injury. Patients with this problem usually display a triad of glomerulonephritis (hematuria), pulmonary hemorrhage (hemoptysis), and antibodies to basement membrane antigens.

Lung or renal biopsies are required to obtain tissue. Serum assays are faster and more reliable methods for diagnosing Goodpasture syndrome.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.

- Fasting: yes

- Blood tube commonly used: red

**PT** If a lung biopsy (see p. 584) or renal biopsy (see p. 776) will be used to collect the specimen, explain these procedures.

### Abnormal findings

#### Positive

Goodpasture syndrome

Autoimmune glomerulonephritis

Lupus nephritis

#### notes

## antiglycan antibodies (Crohn disease prognostic panel, Multiple sclerosis antibody panel)

**Type of test** Blood

### Normal findings

Negative

### Test explanation and related physiology

Antiglycan antibodies are immunologically directed to sugar-containing components on the surface of cells (particularly erythrocytes). Antibodies to glycans can be instigated by bacterial, fungal, and parasite infections. The use of glycan arrays for systematic screening of patients with multiple sclerosis (MS) and inflammatory bowel disease (particularly Crohn disease) has been helpful in differentiating these diseases. Furthermore, these antibodies are used to determine treatment and prognosis.

Using enzyme-linked immunosorbent assays, these antibodies can be identified and quantified. Anti-*Saccharomyces cerevisiae* antibody (ASCA), antilaminaribioside carbohydrate antibody (ALCA), antimannobioside carbohydrate antibody (AMCA), and antichitobioside carbohydrate antibody (ACCA) are used to evaluate Crohn disease and differentiate Crohn colitis from ulcerative colitis. When all are positive, Crohn disease is much more likely than ulcerative colitis.

Other antiglycan antibodies are specific for MS patients.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender, pink, or green

### Abnormal findings

#### ▲ Increased levels

Crohn disease  
Multiple sclerosis

notes

**antiliver/kidney microsomal type 1 antibodies****(Anti-LKM-1 antibodies)**

A

**Type of test** Blood**Normal findings**

≤ 20 units (negative)

20.1-24.9 units (equivocal)

≥ 25 units (positive)

**Test explanation and related physiology**

Autoimmune liver disease (e.g., autoimmune hepatitis and primary biliary cirrhosis) is characterized by the presence of autoantibodies, including smooth muscle antibodies (SMA) (see p. 92), antimitochondrial antibodies (AMA) (see p. 82), and antiliver/kidney microsomal antibodies type 1 (anti-LKM-1). Subtypes of autoimmune hepatitis (AIH) are based on autoantibody reactivity patterns.

Anti-LKM-1 antibodies serve as a serologic marker for AIH type 2 and typically occur in the absence of SMAs and antinuclear antibodies. Patients with AIH type 2 more often tend to be young and female and have a severe form of disease that responds well to immunosuppressive therapy.

Patients with chronic hepatitis resulting from hepatitis C can also have elevated anti-LKM-1 antibodies. The diagnosis of autoimmune liver disease cannot be made on antibody testing alone. In many instances, autoimmune liver disease panel testing is performed.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
  - Fasting: no
  - Blood tube commonly used: serum separator
- PT** Explain the importance of performing the test in the morning.
- Blood may be sent to a reference laboratory. Results are available in about 1 week.

**Abnormal findings****▲ Increased levels**

Autoimmune hepatitis

notes

## antimitochondrial antibody (AMA)

**Type of test** Blood

### Normal findings

No antimitochondrial antibodies (AMAs) at titers  $> 1:5$  or  $< 0.1$  units

### Test explanation and related physiology

The AMA is used primarily to aid in the diagnosis of primary biliary cirrhosis. AMA is an anticytoplasmic antibody directed against a lipoprotein in the mitochondrial membrane. Normally the serum does not contain AMA at a titer greater than a ratio of 1 to 5. AMA appears in 94% of patients with primary biliary cirrhosis. This disease occurs predominantly in young and middle-aged women. It has a slow, progressive course marked by elevated liver enzymes, especially alkaline phosphatase and gamma-glutamyl transpeptidase (see pp. 29 and 435) and positive AMA. Liver biopsy (see p. 571) is usually required to confirm the diagnosis. There are subgroups of AMA. The M-2 subgroup is very specific for primary biliary cirrhosis.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Primary biliary cirrhosis  
 Chronic active hepatitis  
 Systemic lupus erythematosus  
 Syphilis  
 Drug-induced cholestasis  
 Autoimmune hepatitis (e.g., scleroderma, systemic lupus erythematosus)  
 Extrahepatic obstruction  
 Acute infectious hepatitis

notes

## antimyocardial antibody

**Type of test** Blood

### Normal findings

Negative (if positive, serum will be titrated)

### Test explanation and related physiology

This test is used to detect an autoimmune source of myocardial injury and disease. Antimyocardial antibodies (AMAs) may be detected in rheumatic heart disease, cardiomyopathy, postthoracotomy syndrome, and postmyocardial infarction syndromes. This test is used both in the detection of an autoimmune cause for these conditions and for monitoring their response to treatment. Antibodies against heart muscle are also found in 20% to 40% of postcardiac surgery patients and in a smaller number of postmyocardial infarction patients. These antibodies are usually associated with a pericarditis that follows the myocardial injury associated with cardiac surgery or myocardial infarction (Dressler syndrome). AMA has also been detected with cardiomyopathy.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Rheumatic heart disease  
Cardiomyopathy  
Postthoracotomy syndrome  
Postmyocardial infarction  
Rheumatic fever  
Streptococcal infection

notes

## antineutrophil cytoplasmic antibody (ANCA)

**Type of test** Blood**Normal findings**

Components	Reference interval
Antineutrophil cytoplasmic antibody, IgG	<1:20: Not significant
Myeloperoxidase antibody	Negative: $\leq 19$ AU/mL Equivocal: 20-25 AU/mL Positive: $\geq 26$ AU/mL
Serine protease 3 antibody	Negative: $\leq 19$ AU/mL Equivocal: 20-25 AU/mL Positive: $\geq 26$ AU/mL

**Test explanation and related physiology**

ANCAs are directed against cytoplasmic components of neutrophils. This test is used to assist in the diagnosis of granulomatous vascular diseases, such as Wegener granulomatosis (WG). It also is useful in tracking the course of the disease, monitoring the response to therapy, and providing early detection of relapse. WG is a regional systemic vasculitis in which the small arteries of the kidneys, lungs, and upper respiratory tract (nasopharynx) are damaged by a granulomatous inflammation. Diagnosis used to be made by biopsy of clinically affected tissue. Serologic testing now plays a key role in the diagnosis of WG and other systemic vasculitis syndromes.

When ANCAs are detected with indirect immunofluorescence microscopy, two major patterns of staining are present: cytoplasmic ANCA (c-ANCA) and perinuclear ANCA (p-ANCA). Specific immunochemical assays demonstrate that c-ANCA consists mainly of antibodies to proteinase 3 (PR3), and p-ANCA consists of antibodies to myeloperoxidase (MPO). Using the antigen-specific immunochemical assay to characterize ANCA (rather than the pattern of immunofluorescence microscopy) is more specific and more clinically relevant; therefore the terms *proteinase 3-ANCA (PR3-ANCA)* and *myeloperoxidase-ANCA (MPO-ANCA)* are used.

The PR3 autoantigen is highly specific (95%–99%) for WG. When the disease is limited to the respiratory tract, the PR3 is positive in about 65% of patients. Nearly all patients with WG limited to the kidney do not have positive PR3. When WG is inactive, the percentage of positive PR3 drops to about 30%.



The MPO autoantigen is found in 50% of patients with WG centered in the kidney. It also occurs in patients with non-WG glomerulonephritis, such as microscopic polyangiitis (MPA). P-ANCA antibodies can also differentiate various forms of inflammatory bowel disease. (See also antiglycan antibodies, p. 80). P-ANCA antibodies are found in 50% to 70% of patients with ulcerative colitis (UC) but in only 20% of patients with Crohn disease (CD).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: verify with laboratory

### Abnormal findings

#### ▲ Increased levels

Wegener granulomatosis  
Microscopic polyarteritis  
Idiopathic crescentic glomerulonephritis  
Ulcerative colitis  
Primary sclerosing cholangitis  
Autoimmune hepatitis  
Churg-Strauss vasculitis  
Active viral hepatitis  
Crohn disease

### notes

**antinuclear antibody (ANA)****Type of test** Blood**Normal findings**

Negative at 1:40 dilution

**Test explanation and related physiology**

ANA is a group of antinuclear antibodies used to diagnose systemic lupus erythematosus (SLE) and other autoimmune (rheumatic) diseases (Box 1). Some of the antibodies in this group are specific for SLE, and others are specific for other autoimmune diseases. ANA can be tested as a specific antibody or as a group with nonspecific antigens (Box 2). The former is more specific, but testing ANA with less specific antigens may be an excellent preliminary test for those suspected of having autoimmune diseases.

Because almost all patients with SLE develop autoantibodies, a negative ANA test result excludes the diagnosis. Positive results occur in approximately 95% of patients with this disease; however, many other rheumatic diseases (see Box 1) are also associated with ANA.

ANA tests are performed using different assays (indirect immunofluorescence microscopy or ELISA), and results are reported as a titer with a particular type of immunofluorescence pattern (when positive). Whereas low-level titers are considered negative, increased titers are positive and indicate an elevated concentration of antinuclear antibodies.

**BOX 1 Diseases associated with antinuclear antibodies**

- Systemic lupus erythematosus
- Sjögren syndrome
- Scleroderma
- Raynaud disease
- Rheumatoid arthritis
- Dermatomyositis
- Mixed connective tissue disease
- Autoimmune hepatitis
- Autoimmune thyroiditis
- Juvenile rheumatoid arthritis
- Primary biliary cirrhosis
- Polymyositis

**BOX 2 Antinuclear antibodies**


Anti-RNA antibodies  
 Anti-ENA antibodies  
   Antinuclear Smith (SM) antibodies  
   Antinuclear RNP antibodies  
   Antihistidyl antibodies  
 Antichromatin antibodies  
   Antinucleosome antibodies  
   Antihistone antibodies  
 Anti-DNA antibodies  
   Anti-ds-DNA antibodies  
   Anti-ss-DNA antibodies  
 Anti-SS-A (Ro)  
 Anti-SS-B (La)

ANA shows up on indirect immunofluorescence as fluorescent patterns in cells that are fixed to a slide and are evaluated under a UV microscope. Different patterns are associated with a variety of autoimmune disorders. When combined with a more specific subtype of ANA (see [Box 2](#)), the pattern can increase specificity of the ANA subtypes for the various autoimmune diseases ([Figure 3](#)). An example of a positive result might be: “Positive at 1:320 dilution with a homogeneous pattern.”

As the disease becomes less active because of therapy, the ANA titers can be expected to fall ([Table 1](#)). In this text, the more commonly used ANA subtypes are separately discussed.

**TABLE 1 Autoimmune disease and positive ANAs**

Autoimmune disease	Positive antibodies
SLE	ANA, SLE prep, ds-DNA, ss-DNA, anti-DNP, SS-A
Drug-induced SLE	ANA
Sjögren syndrome	RF, ANA, SS-A, SS-B
Scleroderma	ANA, Scl-70, RNA, ds-DNA
Raynaud disease	ACA, Scl-70
Mixed connective tissue disease	ANA, RNP, RF, ss-DNA
Rheumatoid arthritis	RF, ANA, RANA, RAP
Primary biliary cirrhosis	AMA
Thyroiditis	Antimicrosomal, antithyroglobulin
Chronic active hepatitis	ASMA

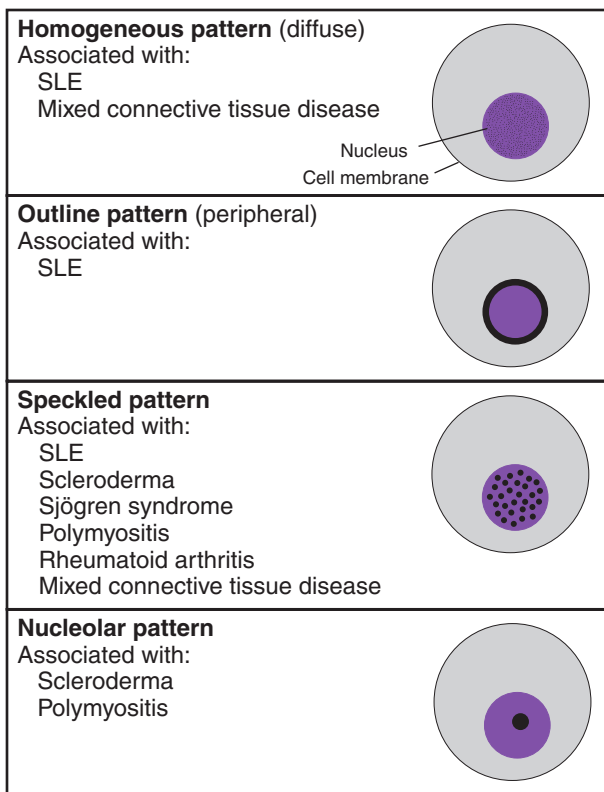

**FIG. 3** Patterns of immunofluorescent staining of antinuclear antibodies and the diseases with which they are associated. *SLE*, Systemic lupus erythematosus.

About 95% of SLE patients have a positive ANA test result. If a patient also has symptoms of SLE (e.g., arthritis, rash, autoimmune thrombocytopenia), then he or she probably has SLE.

### Interfering factors

- ✦ Drugs that may cause a false-positive ANA test include acetazolamide, aminosalicic acid, chlorprothixene, chlorothiazide, griseofulvin, hydralazine, penicillin, phenylbutazone, phenytoin sodium, procainamide, streptomycin, sulfonamides, and tetracyclines.
- ✦ Drugs that may cause a false-negative test include steroids.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

## Abnormal findings

### ▲ Increased levels

Systemic lupus erythematosus  
Rheumatoid arthritis  
Chronic hepatitis  
Periarteritis (polyarteritis) nodosa  
Dermatomyositis  
Scleroderma  
Infectious mononucleosis  
Raynaud disease  
Sjögren syndrome  
Other immune diseases  
Leukemia  
Myasthenia gravis  
Cirrhosis

## notes

**antiparietal cell antibody (APCA)****Type of test** Blood**Normal findings**

Negative

**Test explanation and related physiology**

Parietal cells exist in the proximal stomach and produce hydrochloric acid and intrinsic factor. Intrinsic factor is necessary for the absorption of vitamin B<sub>12</sub> (see p. 969). Antiparietal cell antibodies (APCAs) are found in nearly 90% of patients with pernicious anemia. Nearly 60% of these patients also have antiintrinsic factor antibodies. It is thought that these antibodies contribute to the destruction of the gastric mucosa. APCA is also found in patients with atrophic gastritis, gastric ulcers, and gastric cancer.

APCA is present in other autoimmune-mediated diseases such as thyroiditis, myxedema, juvenile diabetes, Addison disease, and iron-deficiency anemia. Nearly 10% to 15% of the normal population has APCA. As one ages, the incidence of having APCA increases (especially in relatives of patients with pernicious anemia).

APCA can crossreact with other antibodies, especially anticellular and antithyroid antibodies. Titer levels greater than 1:240 are considered positive.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**Abnormal findings****▲ Increased levels**

Pernicious anemia  
Atrophic gastritis  
Hashimoto thyroiditis  
Myxedema  
Insulin-dependent diabetes mellitus  
Addison disease

**notes**

## antiscleroderma antibody (Scl-70 antibody, Scleroderma antibody, RNA polymerase III antibody)

**Type of test** Blood

### Normal findings

Negative

### Test explanation and related physiology

This antibody is diagnostic for scleroderma (progressive systemic sclerosis [PSS]) and is present in 45% of patients with that disease. Scl-70 antibody is an *antinuclear antibody* (see p. 86). PSS is a multisystem disorder characterized by inflammation with subsequent fibrosis of the small blood vessels in skin and visceral organs, including the heart, lungs, kidneys, and gastrointestinal tract. A collagen-like substance is also deposited into the tissue of these organs.

The absence of this antibody does not exclude the diagnosis of PSS. The antibody is rather specific for PSS but is occasionally seen in other autoimmune diseases.

RNA polymerase III antibodies are found in 11% to 23% of patients with PSS. PSS patients who are positive for RNA polymerase III antibodies form a distinct serologic subgroup and usually do not have any of the other antibodies typically found in PSS patients, such as anticentromere (see p. 67) or anti-Scl70. PSS patients with anti-RNA polymerase III have an increased risk of the diffuse cutaneous form of scleroderma, with a high likelihood of skin involvement and hypertensive renal disease.

### Interfering factors

- ✦ Drugs that may cause *increased* levels include aminosalicic acid, isoniazid, methyldopa, penicillin, propylthiouracil, streptomycin, and tetracycline.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### Positive results

Scleroderma

CREST syndrome

notes

## antismooth muscle antibody (ASMA)

**Type of test** Blood

### Normal findings

No antismooth muscle antibodies (ASMAs) at titers > 1:20

### Test explanation and related physiology

The ASMA is used primarily to aid in the diagnosis of autoimmune chronic active hepatitis (CAH), which has also been referred to as *lupoid* CAH. ASMA is an *anticytoplasmic antibody* directed against actin, a cytoskeletal protein. Normally the serum does not contain ASMA at a titer greater than 1:20. ASMA is the most commonly recognized autoantibody in the setting of CAH. It appears in 70% to 80% of patients with CAH. Some types of CAH do not have positive ASMA antibodies.

ASMA is not specific for CAH and can be positive in patients with viral infections, malignancy, multiple sclerosis, primary biliary cirrhosis, and *Mycoplasma* infections. Usually the titer of ASMA is low in these diseases. With CAH, the titer is usually higher than 1 to 160. The titers are not helpful in prognosis nor do they indicate disease response to therapy.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Chronic active hepatitis  
 Mononucleosis hepatitis  
 Primary biliary cirrhosis  
 Viral hepatitis  
 Multiple sclerosis  
 Malignancy  
 Intrinsic asthma

### notes



**antispermatozoal antibody** (Sperm agglutination and inhibition, Sperm antibodies, Antisperm antibodies, Infertility screen)

**Type of test** Fluid analysis; blood

**Normal findings**

< 50% binding

**Test explanation and related physiology**

The antispermatozoal antibody test is an infertility test used to detect the presence of sperm antibodies. Antibodies directed toward sperm antigens can result in diminished fertility. This test is used in the evaluation of an infertile couple usually after a postcoital test result is positive. For fertilization to occur, the sperm head must first attach to the *zona pellucida* of the egg. Sperm antibodies interfere with this binding. Although there is consensus that these antibodies play a role in infertility, the percentage of sperm that must be bound by antibodies before fertility is adversely affected is less clear. IgA antisperm antibodies attached to the sperm tail are associated with poor motility and poor penetration of cervical mucus. IgG antisperm antibodies are associated with blockage of sperm-ovum fusion. Semen and serum may contain sperm antibodies. Semen is the preferred specimen type for men. In cases in which semen production may present difficulties, a serum specimen can be tested instead. Serum is the preferred specimen type in females.

Positives are reported as percentage of sperm with positive bindings, the class of antibody involved (IgG, IgA, and IgM), and the site of binding (head, midpiece, tail, and/or tail tip). Greater than 50% binding is usually required to significantly lower a patient's fertility.

Not only is this test indicated for male infertility studies, but it is also used as a follow-up test when sperm agglutination is noted in the ejaculate. It is also used in men with a history of testicular trauma, biopsy, vasectomy reversal, genital tract infection, or obstructive lesions of the male ductal system. Antisperm antibodies may be found in the blood of men with blocked efferent ducts of the testes (a common cause of low sperm counts or poor sperm mobility) and in 30% to 70% of men who have had a vasectomy. Resorption of sperm from the blocked ducts results in the formation of autoantibodies to sperm as a result of sperm antigens interacting with the immune system. High titers of IgG autoantibodies are often associated with postvasectomy degeneration of the testes, which explains why 50% of males remain infertile after successful repair of a previous vasectomy.

## 94 antispermatozoal antibody

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Sperm specimen

- PT** Inform the man that a semen specimen should be collected after avoiding ejaculation for at least 3 days.
- Give the male patient the proper container for the sperm collection.
- PT** If the specimen is to be collected at home, be certain the patient is told that it must be taken to the laboratory for testing within 2 hours after collection.
- Collect venous blood samples from both the male and the female patient.
- For a vaginal mucus specimen, collect 1 mL of cervical mucus and place it in a plastic vial.
- PT** Instruct the couple on when and how to obtain the test results.

### Abnormal findings

Infertility

Blocked efferent ducts in the testes

Vasectomy

Testicular trauma

notes

## anti-SS-A (RO), anti-SS-B (LA), and anti-SS-C antibodies (Anti-Ro, Anti-La, Sjögren antibodies)

A

**Type of test** Blood

### Normal findings

SS-A (Ro) antibodies, IgG

< 1 U (negative)

≥ 1 U (positive)

SS-B (La) antibodies, IgG

< 1 U (negative)

≥ 1 U (positive)

### Test explanation and related physiology

These three antinuclear antibodies are considered *antiextractable nuclear antigens* (see p. 75) and are used to diagnose Sjögren syndrome. Ro, La, and SS-C antibodies are subtypes of *antinuclear antibodies* (ANAs); they react to nuclear antigens extracted from human B lymphocytes. Ro and La result in a speckled immunofluorescent pattern when seen in the UV microscope (see [Figure 3](#), p. 88). Sjögren syndrome is an immunogenic disease characterized by progressive destruction of the lacrimal and salivary exocrine glands, leading to mucosal and conjunctival dryness. This disease can occur by itself (primary) or in association with other autoimmune diseases such as systemic lupus erythematosus (SLE), rheumatoid arthritis (RA), and scleroderma. In the latter case, it is referred to as secondary Sjögren syndrome.

Anti-SS-A antibodies may be found in approximately 60% to 70% of patients with primary Sjögren syndrome. Anti-SS-B antibodies may be found in approximately half of patients with primary Sjögren syndrome. When anti-SS-A and anti-SS-B antibodies are both positive, Sjögren syndrome can be diagnosed. These antibodies are only occasionally found when secondary Sjögren syndrome is associated with RA. In fact, anti-SS-B is found only in primary Sjögren syndrome. However, anti-SS-C is positive in about 75% of patients with RA or patients with RA and secondary Sjögren syndrome. Therefore these antibodies are also useful in differentiating primary from secondary Sjögren syndrome.

Anti-SS-A can also be found in 25% of patients with SLE. This is particularly useful in ANA-negative cases of SLE because these antibodies are present in the majority of such patients. Anti-SS-B is rarely found in SLE, however. In general, the higher the titer

## 96 anti-SS-A (RO), anti-SS-B (LA), and anti-SS-C antibodies

of anti-SS antibodies, the more likely that Sjögren syndrome exists and the more active the disease is. As Sjögren syndrome becomes less active with therapy, the anti-SS antibody titers can be expected to fall.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### *Positive*

Sjögren syndrome

Rheumatoid arthritis

ANA-negative systemic lupus erythematosus

Neonatal lupus

### notes

**antithrombin activity and antigen assay** (Antithrombin III [AT-III] activity/assay, Functional antithrombin III assay, Heparin cofactor, Immunologic antithrombin III, Serine protease inhibitor)

**Type of test** Blood

**Normal findings**

Antithrombin activity:

Newborn: 35%-40%

Older than 6 months to adult: 80%-130%

Antithrombin antigen assay:

Plasma: > 50% of control value

Serum: 15%-34% lower than plasma value

Immunologic: 17-30 mg/dL

Functional: 80%-120%

Values vary according to laboratory methods.

**Test explanation and related physiology**

AT-III is an  $\alpha_2$  globulin produced in the liver. It inhibits the serine proteases involved in coagulation (II, X, IX, XI, XII). In normal homeostasis, coagulation results from a balance between AT-III and thrombin. A deficiency of AT-III increases coagulation or the tendency toward thrombosis. A hereditary deficiency of AT-III is characterized by a predisposition toward thrombus formation. This is passed on as an autosomal dominant abnormality. Individuals with hereditary AT-III deficiency typically develop thromboembolic events in their early twenties. These thrombotic events are usually venous.

Acquired AT-III deficiency may be seen in patients with cirrhosis, liver failure, advanced carcinoma, nephrotic syndrome, disseminated intravascular coagulation (DIC), protein-losing enteropathies, and acute thrombosis. AT-III is also decreased as much as 30% in pregnant women and women who take estrogens. Antithrombin activity testing is ordered, along with other tests for hypercoagulable disorders (e.g., protein C and protein S, and lupus anticoagulant), when a patient has been experiencing recurrent venous thrombosis.

AT-III provides most of the anticoagulant effect of heparin. Heparin increases antithrombin activity by a thousandfold. Patients who are deficient in AT-III may be heparin resistant and require unusually high doses for an anticoagulation effect. In general, patients respond to heparin if more than 60% of normal AT-III levels exist.

## 98 antithrombin activity and antigen assay

There are two tests for AT-III. The first is a *functional* assay and measures AT-III activity. The second quantifies the AT-III antigen. The antithrombin activity test is performed before the antigen test to evaluate whether the total amount of functional antithrombin activity is normal. Antithrombin activity is the primary (screening) antithrombin assay. If antithrombin activity is normal, AT-III is not the cause of the hypercoagulable state. If antithrombin activity is abnormal, antithrombin antigen should be quantified.

### Interfering factors

- ▼ Drugs that may cause *increased* levels include anabolic steroids, androgens, oral contraceptives (containing progesterone), and sodium warfarin.
- ▼ Drugs that may cause *decreased* levels include fibrinolytic, heparin, L-asparaginase, and oral contraceptives (containing estrogen).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: light blue or red
- Patients receiving heparin therapy may develop a hematoma at the venipuncture site.

### Test results and clinical significance

#### ▲ Increased levels

Kidney transplant  
Acute hepatitis  
Obstructive jaundice  
Vitamin K deficiency

#### ▼ Decreased levels

Disseminated intravascular coagulation (DIC)  
Hypercoagulation states (e.g., deep vein thrombosis)  
Hepatic disorders (especially cirrhosis)  
Nephrotic syndrome  
Protein-wasting diseases (malignancy)  
Hereditary familial deficiency of AT-III

## antithyroglobulin antibody (Thyroid autoantibody, Thyroid antithyroglobulin antibody, Thyroglobulin antibody)

**Type of test** Blood

### Normal findings

< 116 IU/mL

### Test explanation and related physiology

This test is used as a marker for autoimmune thyroiditis and related diseases. Thyroglobulin autoantibodies bind thyroglobulin (Tg), which is a major thyroid-specific protein that plays a crucial role in thyroid hormone synthesis, storage, and release. Tg remains in the thyroid follicles until hormone production is required. Tg is not secreted into the systemic circulation under normal circumstances. However, follicular destruction through inflammation (Hashimoto thyroiditis or chronic lymphocytic thyroiditis and autoimmune hypothyroidism), hemorrhage (nodular goiter), or rapid disordered growth of thyroid tissue (as may be observed in Graves disease or follicular cell-derived thyroid neoplasms) can result in leakage of Tg into the bloodstream. This results in the formation of autoantibodies to Tg in some individuals.

The anti-Tg test is usually performed in conjunction with the antithyroid peroxidase antibody test (see p. 100).

### Interfering factors

- Normal individuals, especially elderly women, may have anti-Tg antibodies.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Hashimoto thyroiditis  
 Rheumatoid arthritis  
 Rheumatoid-collagen disease  
 Pernicious anemia  
 Thyrotoxicosis  
 Hypothyroidism  
 Thyroid carcinoma  
 Myxedema  
 Autoimmune hemolytic anemia

## antithyroid peroxidase antibody (Anti-TPO, TPO-Ab, Antithyroid microsomal antibody, Thyroid autoantibody)

**Type of test** Blood

### Normal findings

Titer < 9 IU/mLco

### Test explanation and related physiology

This test is primarily used in the differential diagnosis of thyroid diseases. Thyroid microsomal antibodies are commonly found in patients with various thyroid diseases. They are present in 70% to 90% of patients with Hashimoto thyroiditis. Microsomal antibodies are produced in response to microsomes escaping from the thyroid epithelial cells surrounding the thyroid follicle. These escaped microsomes then act as antigens and stimulate the production of antibodies. These immune complexes initiate inflammatory and cytotoxic effects on the thyroid follicle.

Although many different thyroid diseases are associated with elevated antimicrosomal antibody levels, the most frequent is chronic thyroiditis (Hashimoto thyroiditis in adults and lymphocytic thyroiditis in children and young adults).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Hashimoto thyroiditis  
 Rheumatoid arthritis  
 Rheumatoid-collagen disease  
 Pernicious anemia  
 Thyrotoxicosis  
 Hypothyroidism  
 Thyroid carcinoma  
 Myxedema

notes



**apolipoproteins**

**Type of test** Blood

**Normal findings****Apo A-I**

Adult/elderly:

Male: 75-160 mg/dL

Female: 80-175 mg/dL

Child:

Newborn:

Male: 41-93 mg/dL

Female: 38-106 mg/dL

6 months-4 years:

Male: 67-167 mg/dL

Female: 60-148 mg/dL

5-17 years: 83-151 mg/dL

**Apo B**

Adult/elderly:

Male: 50-125 mg/dL

Female: 45-120 mg/dL

Child:

Newborn: 11-31 mg/dL

6 months-3 years: 23-75 mg/dL

5-17 years:

Male: 47-139 mg/dL

Female: 41-132 mg/dL

**Apo A-I/Apo B ratio**

Male: 0.85-2.24

Female: 0.76-3.23

**Lipoprotein (a)**

White (5th-95th percentiles):

Male: 2.2-49.4 mg/dL

Female: 2.1-57.3 mg/dL

African American (5th-95th percentiles):

Male: 4.6-71.8 mg/dL

Female: 4.4-75 mg/dL

**Test explanation and related physiology**

This test is used to evaluate the risk of atherogenic heart and peripheral vascular diseases. These proteins are indicators of atherogenic risks much like high-density lipoprotein (HDL),

## 102 apolipoproteins

low-density lipoprotein (LDL), and very-low-density lipoprotein (VLDL) (see p. 565). Apolipoproteins are the protein part of lipoproteins (e.g., HDL, LDL). In general, apolipoproteins play an important role in lipid transport in the lymphatic and circulatory systems. They act as enzyme cofactors in lipoprotein synthesis. Apolipoproteins also act as receptor ligands to improve transport of fat particles in the cell. Apolipoprotein synthesis in the liver is controlled by many factors, including dietary composition, hormones (e.g., insulin, glucagon, thyroxin, estrogens, androgens), alcohol intake, and various drugs (e.g., statins, niacin, fibric acids).

There are several types of apolipoproteins, including apo A-I, apo B, and apo E. *Apolipoprotein A (apo A)* is the major polypeptide component of HDL. Low levels of apo A are associated with increased risk of coronary or peripheral artery disease (CPAD). Elevated levels may protect against CPAD.

*Apo B* is the major polypeptide component of LDL and chylomicrons. Apo B makes cholesterol soluble for deposition in the arterial wall. Forty percent of the protein portion of VLDL is composed of apo B. Familial hypercholesterolemia type B is caused by mutations in the apo B gene.

*Lp(a)* (referred to as *lipoprotein little a*) is a heterogeneous group of lipoproteins consisting of an apo A molecule attached to an apo B molecule. An increased level of Lp(a) may be an independent risk factor for atherosclerosis and is particularly harmful to the endothelium. Serum concentrations of Lp(a) appear to be largely related to genetic factors; diet and statin drugs do not have a major effect on Lp(a) levels. Niacin does lower Lp(a) levels, however.

*Apolipoprotein E (apo E)* is involved in cholesterol transport. Through genotyping, three alleles for apo E have been identified: E2, E3, and E4. Each person gets an allele from each parent. E3/3 is the normal. E2/2 is found rarely and is associated with type III hyperlipidemia. E4/4 or E4/3 is associated with high LDL levels. The apo E4 gene has been proposed as a risk factor for Alzheimer disease (AD). When diminished in the CSF, the risk of AD is increased and cognitive function is decreased (see amyloid beta protein. p. 58). Apo E2 and E4 are associated with increased triglycerides.

*Lp-PLA2* is a lipase enzyme located on the surface of circulating LDL. This protein is atherogenic.

### Interfering factors

#### Apo A-I

- Physical exercise may increase apo A-I levels.

- Smoking may decrease levels.
- Diets high in carbohydrates or polyunsaturated fats may decrease apo A-I levels.
- ▼ Drugs that may *increase* apo A-I levels include carbamazepine, estrogens, ethanol, lovastatin, niacin, oral contraceptives, phenobarbital, pravastatin, and simvastatin.
- ▼ Drugs that may *decrease* apo A-I levels include androgens, beta-blockers, diuretics, and progestins.

### Apo B

- Diets high in saturated fats and cholesterol may increase apo B levels.
- ▼ Drugs that may *increase* apo B levels include androgens, beta-blockers, diuretics, ethanol abuse, and progestins.
- ▼ Drugs that may *decrease* apo B levels include cholestyramine, estrogen (postmenopausal women), lovastatin, neomycin, niacin, simvastatin, and thyroxine.

### Lipoprotein (a)

- ▼ Drugs that may *decrease* Lp(a) include estrogens, neomycin, niacin, and stanozolol.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes (12-14 hours)
- Blood tube commonly used: red
- PT Inform the patient that smoking is prohibited before the test.

### Abnormal findings

#### ▲ Increased apo A-I

Familial  
hyperalphalipoproteinemia  
Pregnancy  
Weight reduction

#### ▼ Decreased apo A-I

Coronary artery disease  
Ischemic coronary disease  
Myocardial infarction  
Familial  
  hyposalphalipoproteinemia  
Fish eye disease  
Uncontrolled diabetes  
  mellitus  
Tangier disease  
Nephrotic syndrome  
Chronic renal failure  
Cholestasis  
Hemodialysis

## 104 apolipoproteins

### ▲ Increased apo B

- Hyperlipoproteinemia (types IIa, IIb, IV, V)
- Nephrotic syndrome
- Pregnancy
- Hemodialysis
- Biliary obstruction
- Coronary artery disease
- Diabetes mellitus
- Hypothyroidism
- Anorexia nervosa
- Renal failure

### ▲ Increased Lp(a)

- Premature coronary artery disease
- Stenosis of cerebral arteries
- Uncontrolled diabetes mellitus
- Severe hypothyroidism
- Familial hypercholesterolemia
- Chronic renal failure
- Estrogen depletion

### Apo E4 gene

- Alzheimer disease

### ▼ Decreased apo B

- Tangier disease
- Hyperthyroidism
- Malnutrition
- Inflammatory joint disease
- Chronic pulmonary disease
- Weight reduction
- Chronic anemia
- Reye syndrome

### ▼ Decreased Lp(a)

- Alcoholism
- Malnutrition
- Chronic hepatocellular disease

notes

**apt test** (Downey test, Qualitative fetal hemoglobin stool test, Stool for swallowed blood)

**Type of test** Stool

**Normal findings**

No fetal blood present.

Maternal blood may be present.

**Test explanation and related physiology**

Blood in the stool of a newborn must be rapidly evaluated. Furthermore, some serious diseases present as rectal bleeding in newborns. Much more commonly, however, newborns may simply be defecating maternal blood that was swallowed during birth or breastfeeding.

The Apt test is performed on the stool specimen to differentiate maternal from fetal blood in the stool. Fetal hemoglobin is resistant to denaturation; adult hemoglobin (hemoglobin A) is not. This test can be performed on stool, a stool-stained diaper, amniotic fluid, or vomitus.

**Procedure and patient care****Before**

**PT** Explain the procedure to the newborn's parents.

- Assess the vital signs of the newborn with possible intestinal bleeding.

**During**

- Obtain an adequate stool or vomitus specimen.
- In the laboratory, 1% NaOH is added to the specimen. Maternal blood turns brown; newborn blood stays red or pink.

**After**

- If maternal blood is present, reassure the parents and examine the mother for nipple erosion or cracking.
- If newborn blood is present, begin close observation and support during further diagnostic procedures.

**Abnormal findings**

Active gastrointestinal bleeding

Necrotizing enterocolitis

notes

## arterial blood gases (ABGs, blood gases)

**Type of test** Blood

### Normal findings

#### pH

Adult/child: 7.35-7.45

Newborn: 7.32-7.49

2 months-2 years: 7.34-7.46

pH (venous): 7.31-7.41

#### Pco<sub>2</sub>

Adult/child: 35-45 mm Hg

Child < 2 years: 26-41 mm Hg

Pco<sub>2</sub> (venous): 40-50 mm Hg

#### HCO<sub>3</sub>

Adult/child: 21-28 mEq/L

Newborn/infant: 16-24 mEq/L

#### Po<sub>2</sub>

Adult/child: 80-100 mm Hg

Newborn: 60-70 mm Hg

Po<sub>2</sub> (venous): 40-50 mm Hg

#### O<sub>2</sub> saturation

Adult/child: 95%-100%

Elderly: 95%

Newborn: 40%-90%

#### O<sub>2</sub> content

Arterial: 15-22 vol %

Venous: 11-16 vol %

#### Base excess

0 ± 2 mEq/L

#### Alveolar to arterial O<sub>2</sub> difference

< 10 mm Hg

### Possible critical values

pH: <7.25, >7.55

PCO<sub>2</sub>: < 20, > 60

HCO<sub>3</sub>: < 15, > 40

PO<sub>2</sub>: < 40

O<sub>2</sub> saturation: 75% or lower

Base excess: ± 3 mEq/L

## Test explanation and related physiology

Measurement of ABGs provides valuable information in assessing and managing a patient's respiratory (ventilation) and metabolic (renal) acid-base and electrolyte homeostasis. It is also used to assess adequacy of oxygenation. ABGs are used to monitor patients on ventilators, to monitor critically ill nonventilator patients, to establish preoperative baseline parameters, and to enlighten electrolyte therapy.

### pH

The pH is inversely proportional to the actual hydrogen ion concentration. Therefore as the hydrogen ion concentration decreases, the pH increases and vice versa. The pH is a measure of alkalinity (pH > 7.4) and acidity (pH < 7.35). In respiratory or metabolic alkalosis, the pH is elevated; in respiratory or metabolic acidosis, the pH is decreased (Table 2).

### Pco<sub>2</sub>

The Pco<sub>2</sub> is a measure of the partial pressure of CO<sub>2</sub> in the blood. Pco<sub>2</sub> is a measurement of ventilation capability. The faster and more deeply one breathes, the more CO<sub>2</sub> is blown off and Pco<sub>2</sub> levels drop. Therefore Pco<sub>2</sub> is referred to as the *respiratory* component in acid-base determination because this value is controlled primarily by the lungs. As the CO<sub>2</sub> level increases, the pH decreases. The CO<sub>2</sub> level and the pH are inversely proportional. The Pco<sub>2</sub> in the blood and cerebrospinal fluid is a major stimulant to the breathing center in the brain. As Pco<sub>2</sub> levels rise, breathing is stimulated. If Pco<sub>2</sub> levels rise too high, breathing cannot keep up with the demand to blow off or ventilate. As Pco<sub>2</sub> levels rise further, the brain is depressed, and ventilation decreases further, causing coma.

The Pco<sub>2</sub> level is elevated in primary respiratory acidosis and is decreased in primary respiratory alkalosis (see Table 2). Because the lungs compensate for primary metabolic acid-base derangements, Pco<sub>2</sub> levels are affected by metabolic disturbances as well. In metabolic acidosis, the lungs attempt to compensate by blowing off CO<sub>2</sub> to raise pH. In metabolic alkalosis, the lungs attempt to compensate by retaining CO<sub>2</sub> to lower pH (Table 3).

### Bicarbonate (HCO<sub>3</sub>) or CO<sub>2</sub> content

Most of the CO<sub>2</sub> content in the blood is HCO<sub>3</sub>. The bicarbonate ion is a measure of the *metabolic (renal/kidney)* component of the acid-base equilibrium. It is regulated by the kidneys. This ion can be measured directly by the bicarbonate value or indirectly by the CO<sub>2</sub> content (see p. 197). As the HCO<sub>3</sub> level

**TABLE 2** Normal values for arterial blood gases and abnormal values in uncompensated acid-base disturbances

Acid-base disturbances	pH	Pco <sub>2</sub> (mm Hg)	HCO <sub>3</sub> (mEq/L)	Common causes
None (normal values)	7.35–7.45	35–45	22–26	
Respiratory acidosis	↓	↑	Normal	Respiratory depression (drugs, central nervous system trauma) Pulmonary disease (pneumonia, chronic obstructive pulmonary disease, respiratory underventilation)
Respiratory alkalosis	↑	↓	Normal	Hyperventilation (emotions, pain, respirator overventilation)
Metabolic acidosis	↓	Normal	↓	Diabetes, shock, renal failure, intestinal fistula
Metabolic alkalosis	↑	Normal	↑	Sodium bicarbonate overdose, prolonged vomiting, nasogastric drainage



**TABLE 3** Acid-base disturbances and compensatory mechanisms

Acid-base disturbance	Mode of compensation
Respiratory acidosis	Kidneys will retain increased amounts of $\text{HCO}_3^-$ to increase pH.
Respiratory alkalosis	Kidneys will excrete increased amounts of $\text{HCO}_3^-$ to lower pH.
Metabolic acidosis	Lungs blow off $\text{CO}_2$ to raise pH.
Metabolic alkalosis	Lungs retain $\text{CO}_2$ to lower pH.

increases, the pH also increases; therefore the relationship of bicarbonate to pH is directly proportional.  $\text{HCO}_3^-$  is elevated in metabolic alkalosis and decreased in metabolic acidosis (see Table 3). The kidneys also are used to compensate for primary respiratory acid-base derangements. For example, in respiratory acidosis, the kidneys attempt to compensate by resorbing increased amounts of  $\text{HCO}_3^-$ . In respiratory alkalosis, the kidneys excrete  $\text{HCO}_3^-$  in increased amounts to lower pH (see Table 3.)

### $\text{Po}_2$

This is an indirect measure of the oxygen content of arterial blood.  $\text{Po}_2$  is a measure of the tension (pressure) of oxygen dissolved in the plasma. This pressure determines the force of  $\text{O}_2$  to diffuse across the pulmonary alveoli membrane. The  $\text{Po}_2$  level is decreased in patients who:

- Are unable to oxygenate the arterial blood because of  $\text{O}_2$  diffusion difficulties (e.g., pneumonia)
- Have premature mixing of venous blood with arterial blood (e.g., in congenital heart disease)
- Have underventilated and overperfused pulmonary alveoli (Pickwickian syndrome or patients with significant atelectasis)

### $\text{O}_2$ saturation

Oxygen saturation is an indication of the percentage of hemoglobin saturated with  $\text{O}_2$ . When 92% to 100% of the hemoglobin carries  $\text{O}_2$ , the tissues are adequately provided with  $\text{O}_2$ , assuming normal  $\text{O}_2$  dissociation. As the  $\text{Po}_2$  level decreases, the percentage of hemoglobin saturation also decreases. When the  $\text{Po}_2$  level drops to lower than 60 mm Hg, small decreases in the  $\text{Po}_2$  level cause large decreases in the percentage of hemoglobin saturated with  $\text{O}_2$ . At  $\text{O}_2$  saturation levels of 70% or lower, the tissues are unable to extract enough  $\text{O}_2$  to carry out their vital functions.

## 110 arterial blood gases

O<sub>2</sub> saturation is calculated by the blood gas machine using the following formula:

$$\text{Percentage of O}_2 \text{ saturation} = 100 \times \frac{\text{Volume of O}_2 \text{ content Hgb}}{\text{Volume of O}_2 \text{ Hgb capacity}}$$

*Pulse oximetry* is a noninvasive method of determining O<sub>2</sub> saturation (p. 658).

### O<sub>2</sub> content

This is a calculated number that represents the amount of O<sub>2</sub> in the blood. The formula for calculation is:

$$\text{O}_2 \text{ content} = \text{O}_2 \text{ saturation} \times \text{Hgb} \times 1.34 + \text{Po}_2 \times 0.003$$

Nearly all O<sub>2</sub> in the blood is bound to hemoglobin. O<sub>2</sub> content decreases with the same diseases that diminish PO<sub>2</sub>.

### Base excess or deficit

This number is calculated by the blood gas machine by using the pH, PCO<sub>2</sub>, and hematocrit. It represents the amount of buffering anions in the blood. HCO<sub>3</sub> is the largest of these. Others include hemoglobin, proteins, and phosphates. Base excess is a way to take all these anions into account when determining acid-base treatment based on the *metabolic* component. Negative base excess (deficit) indicates a metabolic acidosis (e.g., lactic acidosis). A positive base excess indicates metabolic alkalosis or compensation to prolonged respiratory acidosis.

### Alveolar (A) to arterial (a) O<sub>2</sub> difference (A-a gradient)

This is a calculated number that indicates the difference between alveolar (A) O<sub>2</sub> and arterial (a) O<sub>2</sub>. The normal value is less than 10 mm Hg (torr). If the A-a gradient is abnormally high, there is either a problem in diffusing O<sub>2</sub> across the alveolar membrane (thickened edematous alveoli) or unoxygenated blood is mixing with the oxygenated blood.

### Contraindications

- Patients with a negative Allen test result
- Patients with arteriovenous fistula proximal to the access site
- Patients with severe coagulopathy

### Potential complications

- Occlusion of the artery used for access
- Penetration of other important structures anatomically juxtaposed to the artery (e.g., nerve)

## Interfering factors

- O<sub>2</sub> saturation can be falsely increased with the inhalation of carbon monoxide.
- ✚ Respiration can be inhibited by sedatives or narcotics.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that an arterial puncture is associated with more discomfort than a venous puncture.
- Notify the laboratory before drawing ABGs so that the equipment can be calibrated before the blood sample arrives.
- Perform the *Allen test* to assess collateral circulation.
- To perform the Allen test, make the patient's hand blanch by obliterating both the radial and the ulnar pulses, and then release the pressure over the ulnar artery only. If flow through the ulnar artery is good, flushing will be seen immediately. The Allen test is then positive, and the radial artery can be used for puncture.
- If the Allen test is negative (no flushing), repeat it on the other arm.
- If both arms give a negative result, choose another artery for puncture.

### During

- Note that arterial blood can be obtained from any area of the body in which strong pulses are palpable, usually from the radial, brachial, or femoral artery.
- Cleanse the arterial site.
- Use a small gauge needle to collect the arterial blood in an air-free heparinized syringe.
- After drawing blood, remove the needle and apply pressure to the arterial site for 3 to 5 minutes.
- Expel any air bubbles in the syringe.
- Cap the syringe and gently rotate to mix the blood and heparin.
- Note that an arterial puncture is performed by laboratory technicians, respiratory-inhalation therapists, nurses, or physicians in approximately 10 minutes.

### After

- Place the arterial blood on ice and immediately take it to the chemistry laboratory for analysis.
- Apply pressure or a pressure dressing to the arterial puncture site for 3 to 5 minutes to avoid hematoma formation.
- Assess the puncture site for bleeding. Remember that an artery rather than a vein has been stuck.

### Abnormal findings

▲ **Increased pH (alkalosis)**

*Metabolic alkalosis*

- Hypokalemia
- Hypochloremia
- Chronic and high-volume gastric suction
- Chronic vomiting
- Aldosteronism
- Mercurial diuretics

*Respiratory alkalosis*

- Chronic heart failure
- Cystic fibrosis
- Carbon monoxide poisoning
- Pulmonary emboli
- Shock
- Acute and severe pulmonary diseases
- Anxiety neuroses
- Pain
- Pregnancy

▲ **Increased  $\text{PCO}_2$**

- Chronic obstructive pulmonary disease (COPD) (bronchitis, emphysema)
- Oversedation
- Head trauma
- Overoxygenation in a patient with COPD
- Pickwickian syndrome

▲ **Increased  $\text{HCO}_3$**

- Chronic vomiting
- Chronic and high-volume gastric suction
- Aldosteronism
- Use of mercurial diuretics
- Chronic obstructive pulmonary disease

▼ **Decreased pH (acidosis)**

*Metabolic acidosis*

- Ketoacidosis
- Lactic acidosis
- Severe diarrhea
- Renal failure

*Respiratory acidosis*

- Respiratory failure

▼ **Decreased  $\text{PCO}_2$**

- Hypoxemia
- Pulmonary emboli
- Anxiety
- Pain
- Pregnancy

▼ **Decreased  $\text{HCO}_3$**

- Chronic and severe diarrhea
- Chronic use of loop diuretics
- Starvation
- Diabetic ketoacidosis
- Acute renal failure

- ▲ **Increased  $\text{PO}_2$ ,  
increased  $\text{O}_2$  content**  
Polycythemia  
Increased inspired  $\text{O}_2$   
Hyperventilation

- ▼ **Decreased  $\text{PO}_2$ ,  
decreased  $\text{O}_2$  content**  
Anemias  
Mucus plug  
Bronchospasm  
Atelectasis  
Pneumothorax  
Pulmonary edema  
Adult respiratory distress  
  syndrome  
Restrictive lung disease  
Atrial or ventricular cardiac  
  septal defects  
Emboli  
Inadequate  $\text{O}_2$  in inspired  
  air (suffocation)  
Severe hypoventilation (e.g.,  
  oversedation, neurologic  
  somnia)

notes

**arteriography** (Angiography)

**Type of test** X-ray with contrast dye

**Normal findings**

Normal arterial vasculature

**Test explanation and related physiology**

With the injection of radiopaque contrast material into arteries, blood vessels can be visualized to determine arterial anatomy or vascular disease. Blood flow dynamics, arterial occlusive disease, or vascular anomalies are easily seen. With the use of *digital subtraction angiography (DSA)*, bony structures can be obliterated from the picture. Coronary arteriography is described under cardiac catheterization (see p. 201).

*Renal angiography* permits evaluation of renal artery blood flow dynamics. Arteriosclerotic narrowing (stenosis) of the renal artery is best demonstrated with this study. The angiographic location of the stenotic area is helpful if considering surgical repair of stenting.

*Lower extremity arteriography* allows for accurate identification and location of occlusions within the abdominal aorta and lower extremity arteries. Total or near-total occlusion of the flow of dye is seen in arteriosclerotic vascular occlusive disease. Emboli are seen as total occlusions of the artery. Likewise, arterial traumas, such as lacerations or intimal tears (laceration of the inner arterial lining), appear as total or near-total obstruction of the flow of dye.

Arterial vascular balloon dilation and stenting can be performed if a short-segment arterial stenosis is identified.

**Contraindications**

The following list represents relative contraindications. If the information/therapy is necessary to obtain through arteriography, appropriate steps can be taken to reduce risks in these patients. As in all diagnostic testing, the risks must be weighed against the benefits.

- Patients with allergies to shellfish or iodinated dye
- Patients who are pregnant, unless the benefits outweigh the risks
- Patients with renal disorders, because iodinated contrast is nephrotoxic
- Patients with a bleeding propensity

**Potential complications**

- Allergic reaction to iodinated dye
- Hemorrhage from the arterial puncture site

- Arterial embolism from dislodgment of an arteriosclerotic plaque
- Soft tissue infection around the puncture site
- Renal failure, especially in elderly patients who are chronically dehydrated or have a mild degree of renal failure
- Dissection of the intimal lining of the artery causing complete or partial arterial occlusion
- Pseudoaneurysm development as a result of failure of the puncture site to seal
- Lactic acidosis may occur in patients who are taking metformin. The metformin should not be taken the day of the test to prevent this complication.

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient. Allay any fears and allow the patient to verbalize concerns. See p. xviii for radiation exposure and risks.

- Obtain written and informed consent for this procedure.

**PT** Inform the patient that a warm flush may be felt when the dye is injected.

- Assess the possibility of allergies to iodinated dye.
- Determine whether the patient has been taking anticoagulants.
- Keep the patient NPO for 2 to 8 hours before testing.
- Mark the site of the patient's distal peripheral pulses with a pen before arterial catheterization. This will permit assessment of the peripheral pulses after the procedure.
- If the patient does not have peripheral pulses before arteriography, document that fact so that arterial occlusion will not be suspected on the postangiogram assessment.
- Ensure that the appropriate renal function studies are normal.
- If the patient has diminished renal function, provide IV hydration to minimize further renal damage.

**PT** Instruct the patient to void before the study because the iodinated dye can act as an osmotic diuretic.

**PT** Inform the patient that bladder distention may cause some discomfort during the study.

### During

- Note the following preprocedure steps:
  1. The patient may be sedated before being taken to the angiography room, which is usually within the radiology department.
  2. The patient is placed on the x-ray table in the supine position.
  3. The area around the access artery is shaved and prepared as per protocol and draped in a sterile manner.

## 116 arteriography

4. A catheter is threaded under fluoroscopic visualization to the desired position. Through the catheter, iodinated contrast material is injected.
  5. Cinefluoroscopy is used to visualize the injection in real time.
- Note that this procedure is usually performed by an angiographer (radiologist) in approximately 1 hour.
- PT** During the dye injection, remind the patient that an intense, burning flush may be felt throughout the body but lasts only a few seconds.
- PT** Tell the patient that the most significant discomfort is the puncture necessary for arterial access.
- PT** Remind the patient of the discomfort of lying on a hard x-ray table for a long period.

### After

- After x-ray studies are completed, a pressure dressing is applied to the puncture site.
  - Monitor the patient's vital signs for indications of hemorrhage.
  - Assess the distal arterial pulse in the extremity used for vascular access and compare it with the preprocedure baseline.
  - Observe the arterial puncture site frequently for signs of bleeding or hematoma.
  - Keep the patient on bed rest for up to 8 hours after the procedure to allow for complete sealing of the arterial puncture site. If a vascular closure product is used, the patient may ambulate within 2 hours.
  - Note and compare the color and temperature of the extremity with that of the uninvolved extremity.
- PT** Instruct the patient to drink fluids to prevent dehydration caused by the diuretic action of the dye.
- Evaluate the patient for delayed allergic reaction to the dye.
- PT** Instruct the patient to report any signs of numbness, tingling, pain, or loss of function in the involved extremity.

### Abnormal findings

#### Arteriography of the peripheral vascular system

Arteriosclerotic occlusion

Embolus occlusion

Primary arterial diseases (e.g., fibromuscular dysplasia, Buerger disease)

Aneurysm



## **Kidney arteriography**

Atherosclerotic narrowing of the renal artery

Fibrodysplasia of the renal artery

Renal vascular causes of hypertension

**notes**

## arthrocentesis with synovial fluid analysis

**Type of test** Fluid analysis

### Normal findings

Appearance	Clear, straw colored, no blood
RBC	None
WBC	0-150/mm <sup>3</sup>
WBC differential	
Neutrophils	7%
Lymphocytes	24%
Monocytes	48%
Macrophages	10%
Glucose	Equal to fasting blood glucose
Protein	1-3 dL
LDH	< 25 mg/dL
Uric acid	6-8 mg/dL
Gram stain	Negative

### Test explanation and related physiology

Arthrocentesis is performed to establish the diagnosis of joint infection, arthritis, crystal-induced arthritis (gout and pseudogout), synovitis, or neoplasms involving the joint. This procedure is also used to identify the cause of joint inflammation or effusion and to inject antiinflammatory medications (usually corticosteroids) into a joint space.

Arthrocentesis is performed by inserting a sterile needle into the joint space of the involved joint to obtain synovial fluid for analysis. Synovial fluid is a liquid found in small amounts within the joints. Aspiration (withdrawal of the fluid) may be performed on any major joint such as the knee, shoulder, hip, elbow, wrist, or ankle.

The fluid sample is examined. Normal joint fluid is clear, straw colored, and quite viscous. Viscosity is reduced in patients with inflammatory arthritis. Viscosity can be roughly estimated by forcing some synovial fluid from a syringe. Fluid of normal viscosity forms a “string” more than 5 cm long (*string sign*); fluid of low viscosity as seen in inflammation drips in a manner similar to water. A Gram stain and culture of the fluid is usually performed.

The *mucin clot test* correlates with the viscosity and is an estimation of hyaluronic acid–protein complex integrity. This test is performed by adding acetic acid to joint fluid. The formation of a tight, ropy clot indicates qualitatively good mucin and

the presence of adequate molecules of intact hyaluronic acid. Hyaluronic acid can be directly quantified. The mucin clot is poor in quality and quantity in the presence of an inflammatory joint disease such as rheumatoid arthritis (RA). By itself, synovial fluid should not spontaneously form a fibrin clot (clot without the addition of acetic acid) because normal joint fluid does not contain fibrinogen. If, however, bleeding into the joint (from trauma or injury) has occurred, the synovial fluid will clot.

The synovial fluid glucose value is usually within 10 mL/dL of the fasting serum glucose value. For proper interpretation, the synovial fluid glucose and serum glucose samples should be drawn simultaneously after the patient has fasted for 6 hours. The synovial fluid glucose level falls with increasing severity of inflammation. Although lowest in septic arthritis, a low synovial glucose level also may be seen in patients with rheumatoid arthritis. The synovial fluid is also tested for protein, uric acid, and lactate levels. Increased uric acid levels indicate gout. Increased protein and lactate levels indicate bacterial infection or inflammation.

Cell counts are also performed on the synovial fluid. Normally the joint fluid contains less than 150 WBCs/mm<sup>3</sup> and 2000 RBCs/mL. An increased WBC count with a high percentage of neutrophils (over 75%) supports the diagnosis of acute bacterial infectious arthritis. Leukocytes can also occur in other conditions such as acute gouty arthritis and rheumatoid arthritis. The differential WBC count, however, will indicate monocytosis or lymphocytosis with these later-mentioned diseases.

Bacterial and fungal cultures are usually requested and performed when infection is suspected. The administration of antibiotics before arthrocentesis may diminish growth of bacteria from synovial fluid cultures and confound results. Smears for acid-fast stains for tubercle bacilli are also performed on the synovial fluid. Synovial fluid is also examined under polarized light for the presence of crystals, which permits differential diagnosis between gout and pseudogout.

The synovial fluid is also analyzed for complement levels (see p. 266). Complement levels are decreased in patients with systemic lupus erythematosus, rheumatoid arthritis, or other immunologic arthritis. These decreased joint complement levels are caused by consumption of the complement induced by the antigen-antibody immune complexes within the joint cavity.

One of the most important tests routinely performed on synovial fluid is the microscopic examination for crystals. For example, urate crystals indicate gouty arthritis. Calcium pyrophosphate

## 120 arthrocentesis with synovial fluid analysis

crystals are found in pseudogout. Cholesterol crystals occur in rheumatoid arthritis.

### Contraindications

- Patients with skin or wound infections near the needle puncture because of the risk of sepsis

### Potential complications

- Joint infection
- Hemorrhage in the joint area

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Obtain an informed consent if indicated.
- The physician may request that the patient be kept NPO after midnight on the day of the test.

#### During

- Have the patient lie on his or her back with the joint fully extended.
- Note the following procedural steps:
  1. The skin is locally anesthetized to minimize pain.
  2. The area is aseptically cleansed, and a needle is inserted through the skin and into the joint space.
  3. Fluid is obtained for analysis. If a corticosteroid or other medications (e.g., antibiotics) are to be administered, a syringe containing the steroid preparation is attached to the needle, and the drug is injected.
  4. The needle is removed, and a wrap dressing may be applied to the site.
  5. Sometimes a peripheral venous blood sample is taken to compare chemical tests on the blood with chemical studies on the synovial fluid.
- Note that a physician performs this procedure in an office or at the patient's bedside in approximately 20 minutes.

**PT** Tell the patient that the only discomfort associated with this test is the injection of the local anesthetic.

- Be aware that joint-space pain may worsen after fluid aspiration, especially in patients with acute arthritis.

#### After

**PT** Assess the joint for any pain, fever, or swelling. Teach the patient to look for signs of infection at home.

**PT** Apply ice to decrease pain and swelling and instruct the patient to continue this at home.

- PT** Tell the patient to avoid strenuous use of the joint for the next several days.
- PT** Teach the patient to walk on crutches if indicated.
- PT** Instruct the patient to look for signs of bleeding into the joint (significant swelling, increasing pain, or joint weakness).
- PT** Instruct the patient not to drive until it is approved by the physician.

### **Abnormal findings**

Infection  
Osteoarthritis  
Synovitis  
Neoplasm  
Joint effusion  
Septic arthritis  
Systemic lupus erythematosus  
Rheumatoid arthritis  
Gout  
Pseudogout

### **notes**

## arthroscopy

**Type of test** Endoscopy


### Normal findings

Normal ligaments, menisci, and articular surfaces of the joint

### Test explanation and related physiology

Arthroscopy is an endoscopic procedure that allows direct visual examination of a joint space with a specially designed endoscope (Figure 4). Although this technique can visualize many joints of the body, it is most often used to evaluate the knee for meniscus cartilage or ligament injury. It is also used in the differential diagnosis of acute and chronic disorders of the knee (e.g., arthritic inflammation vs. injury). The shoulder is commonly visualized for injury to the rotator cuff.

Physicians can perform corrective surgery through the endoscope. Diagnostic arthroscopy is performed less often because of the availability and accuracy of MRI of the joints (see p. 599).

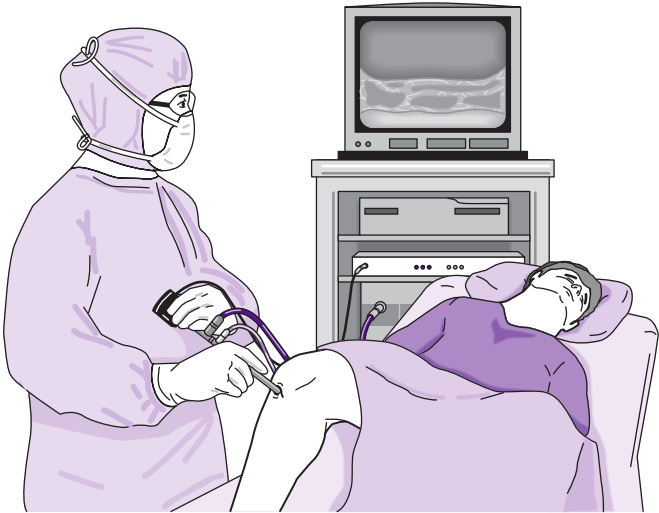

**FIG. 4** Arthroscopy. The arthroscope is placed in the joint space of the knee. Video arthroscopy requires a water source to distend the joint space, a light source to see the contents of the joint, and a television monitor to project the image. Other trocars are used for access of the joint space for other operative instruments.

Visual findings may be recorded by attaching a video camera to the arthroscope. Other joints that can be evaluated by the arthroscope include the tarsal, ankle, knee, hip, carpal, wrist, and temporomandibular joints. Synovial fluid can be obtained for fluid analysis through arthroscopy. See arthrocentesis, p. 118.

### Contraindications

- Patients with ankylosis
- Patients with local skin infections near the joint
- Patients who have recently had an arthrogram

### Potential complications

- Infection
- Hemarthrosis
- Swelling
- Thrombophlebitis
- Joint injury
- Synovial rupture

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Ensure that the physician has obtained written consent for this procedure.
- Follow the routine preoperative procedure of the institution.
- Keep the patient NPO after midnight on the day of the test.

**PT** Instruct the patient who will use crutches after the procedure regarding the appropriate crutch gait. The patient should use crutches after arthroscopy until able to walk without limping.

- Shave the hair in the area 6 inches above and below the joint before the test (if ordered).

#### During

- The patient is placed supine on an operating room table.
- Note the following procedural steps:
  1. Local or general anesthesia is used.
  2. The leg may be scrubbed and wrapped with an elastic bandage from the toes to the lower thigh to drain as much blood from the leg as possible.
  3. A tourniquet is placed on the patient's leg. If the tourniquet is not used, a fluid solution may be instilled into the patient's knee immediately before insertion of the arthroscope to distend the knee and help reduce bleeding.
  4. The arthroscope (a lighted instrument) is inserted into the joint space through a small incision to visualize the inside of the knee joint. Although the entire joint can be

## 124 arthroscopy

viewed from one puncture site, additional punctures may be required for any reparative work.

5. Before removal of the arthroscope, the joint is irrigated. After a few sutures are placed into the skin, a wrapped dressing is applied to the incision site.
- Note that this procedure is performed in the operating room by an orthopedic surgeon in approximately 60 minutes.

### After

- PT** Instruct the patient to elevate the knee when sitting and to avoid overbending the knee, so that swelling is minimized.
- PT** Inform the patient that he or she can usually walk with the assistance of crutches; however, this depends on the extent of the procedure and the physician's protocol.
- PT** Tell the patient to minimize use of the joint for several days.
  - Examine the incision site for bleeding.
- PT** Educate the patient to look for signs of bleeding into the joint (significant swelling, increasing pain, or joint weakness).
  - Apply ice to reduce pain and swelling and instruct the patient to continue this at home.
- PT** Educate the patient to look for signs of phlebitis. This is not uncommon in a person immobilized by joint pain. The involved leg may become swollen, painful, and edematous.
- PT** Instruct the patient not to drive until it is approved by the physician.
- PT** Inform the patient that the sutures will be removed in approximately 7 to 10 days.

### Abnormal findings

Torn cartilage  
Torn ligament  
Patellar disease  
Patellar fracture  
Chondromalacia  
Osteochondritis dissecans  
Cyst (e.g., Baker)  
Synovitis  
Rheumatoid arthritis  
Degenerative arthritis  
Meniscal disease  
Osteochondromatosis  
Trapped synovium

### notes



**aspartate aminotransferase** (AST; Formerly called serum glutamic-oxaloacetic transaminase [SGOT])

**Type of test** Blood

**Normal findings**

Adult: 0-35 units/L or 0-0.58  $\mu$ Kat/L (SI units); females tend to have slightly lower values than males

Elderly: values slightly higher than adult values

Children:

0-5 days: 35-140 units/L

< 3 years: 15-60 units/L

3-6 years: 15-50 units/L

6-12 years: 10-50 units/L

12-18 years: 10-40 units/L

**Test explanation and related physiology**

Because AST exists within the liver cells, diseases that affect the hepatocytes cause elevated levels of this enzyme. This test is used in the evaluation of suspected hepatocellular diseases. The amount of elevation depends on the time after the injury that the blood is drawn. AST is cleared from the blood in a few days. Serum AST levels become elevated 8 hours after cell injury, peak at 24 to 36 hours, and return to normal in 3 to 7 days. If the cellular injury is chronic, levels will be persistently elevated.

In acute hepatitis, AST levels can rise to 20 times the normal value. In acute extrahepatic obstruction (e.g., gallstones), AST levels quickly rise to 10 times the normal value and fall swiftly. In cirrhotic patients, the level of AST depends on the amount of active inflammation.

Serum AST levels are often compared with alanine aminotransferase (ALT, see p. 125) levels. The AST-to-ALT ratio is usually greater than 1.0 in patients with alcoholic cirrhosis, liver congestion, or metastatic tumor of the liver. Ratios less than 1.0 may be seen in patients with acute hepatitis, viral hepatitis, or infectious mononucleosis. The ratio is less accurate if AST levels exceed 10 times the normal value.

Patients with acute pancreatitis, acute renal diseases, musculoskeletal diseases, or trauma may have a transient rise in serum AST. Patients with red blood cell abnormalities, such as acute hemolytic anemia and severe burns, also can have elevations of this enzyme.

## Interfering factors

- Exercise may cause increased levels.
- Pyridoxine deficiency (beriberi or pregnancy), severe long-standing liver disease, uremia, or diabetic ketoacidosis may cause decreased levels.
- ▼ Drugs that may cause *increased* levels include antihypertensives, cholinergic agents, coumarin-type anticoagulants, digitalis preparations, erythromycin, hepatotoxic medications, isoniazid, methyldopa, opiates, oral contraceptives, salicylates, statins, and verapamil.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- If possible, avoid giving the patient any IM injection because increased enzyme levels may result.
- Record the time and date of any IM injection given.
- Record the exact time and date when the blood test is performed. This aids in the interpretation of the temporal pattern of enzyme elevations.

## Abnormal findings

### ▲ Increased levels

#### *Liver diseases*

- Hepatitis
- Hepatic cirrhosis
- Drug-induced liver injury
- Hepatic metastasis
- Hepatic necrosis (initial stages only)
- Hepatic surgery
- Infectious mononucleosis with hepatitis
- Hepatic infiltrative process (e.g., tumor)

#### *Skeletal muscle diseases*

- Skeletal muscle trauma
- Recent noncardiac surgery
- Multiple traumas
- Severe, deep burns

### ▼ Decreased levels

- Acute renal disease
- Beriberi
- Diabetic ketoacidosis
- Pregnancy
- Chronic renal dialysis

**▲ Increased levels***Skeletal muscle diseases**(cont.)*Progressive muscular  
dystrophy

Recent convulsions

Heat stroke

Primary muscle diseases

(e.g., myopathy, myositis)

*Other diseases*

Acute hemolytic anemia

Acute pancreatitis

**notes**

**barium enema** (BE, Lower GI series)

**Type of test** X-ray with contrast dye

**Normal findings**

Normal filling, contour, patency, and positioning of barium in the colon

Normal filling of the appendix and terminal ileum

**Test explanation and related physiology**

The BE study consists of a series of x-rays with contrast to visualize the colon. It is used to demonstrate the presence and location of polyps, tumors, and diverticula. Anatomic abnormalities (e.g., malrotation) also can be detected. Therapeutically, the BE may be used to reduce nonstrangulated ileocolic intussusception in children.

By reflux of barium in the terminal ileum, Crohn disease (regional enteritis) can be identified. Inflammatory bowel disease involving the colon can be detected with a BE. Fistulas involving the colon can be demonstrated by a BE.

In many instances, air is insufflated into the colon after the instillation of barium. This provides an air contrast to the barium. With air contrast, the colonic mucosa can be much more accurately visualized. This is called an *air-contrast BE* and is more accurate than single contrast barium enema.

**Contraindications**

- Patients suspected of a perforation of the colon
- In these patients, diatrizoate (Gastrografin), a water-soluble contrast medium, is used.
- Patients with megacolon

**Potential complications**

- Colonic perforation, especially when the colon is weakened by inflammation, tumor, or infection
- Barium fecal impaction

**Interfering factors**

- Barium within the abdomen from previous barium tests
- Significant residual stool within the colon may be confused with polyps.

**Procedure and patient care****Before**

**PT** Explain the procedure to the patient. Encourage the patient to verbalize questions and fears. See p. xviii for radiation exposure and risks.

- Assist the patient with the bowel preparation, which varies among institutions. In elderly patients, this preparation can be exhausting and may even cause severe dehydration. A typical preparation for most adults would include the following actions:

#### *Day before examination*

- Give the patient clear liquids for lunch and supper (no dairy products).
- PT** Instruct the patient to drink one glass of water or clear fluid every hour for 8 to 10 hours.
- Administer a cathartic (10 oz of magnesium citrate) or X-Prep (extract of senna fruit) at 2 PM. In children, lesser volumes may be used.
  - Administer three 5-mg bisacodyl (Dulcolax) tablets at 7 PM.
  - A pediatric Fleet enema the night before testing and repeated 3 hours before testing may be adequate prep for an infant.
  - Keep the patient NPO (nothing by mouth) after midnight the day of the test.

#### *Day of examination*

- Keep the patient NPO.
  - Administer a bisacodyl suppository at 6 AM and/or a cleansing enema.
  - Note that pediatric patients will have individualized bowel preparations.
  - Note that special preparations will be ordered for patients with an ileostomy or colostomy.
  - Determine whether the bowel is adequately cleansed. When the fecal return is similar to clear water, preparation is adequate; if large, solid fecal waste is still being evacuated, preparation is inadequate. Notify the radiologist, who may want to extend the bowel preparation.
- PT** Suggest that the patient take reading material to the x-ray department to occupy the time while expelling the barium.

#### **During**

- Note the following procedural steps:
  1. The test begins with placement of a rectal balloon catheter.
  2. The balloon on the catheter is inflated tightly against the anal sphincter to hold the barium within the colon.
  3. The patient is asked to roll into the lateral, supine, and prone positions.
  4. The barium is dripped into the rectum by gravity.
  5. The barium flow is monitored fluoroscopically.
  6. The colon is thoroughly examined as the barium flow progresses through the large colon and into the terminal ileum.

## 130 barium enema

7. The barium is drained out.
  8. If an air-contrast BE has been ordered, air is insufflated into the large bowel.
  9. The patient is asked to expel the barium, and a postevacuation x-ray image is taken.
- Note that this test is usually performed in the radiology department by a radiologist in approximately 45 minutes.
- PT** Inform the patient that abdominal bloating and rectal pressure will occur during instillation of barium.

### After

- Ensure that the patient defecates as much barium as possible.
- PT** Suggest the use of soothing ointments on the anal area to minimize any anorectal pain that may result from the test preparation.
- PT** Encourage ingestion of fluids to avoid dehydration caused by the cathartics.
- PT** Inform the patient that bowel movements will be white. When all the barium has been expelled, the stool will return to a normal color.
- Note that laxatives may be ordered to facilitate evacuation of barium.

### Abnormal findings

Malignant tumor

Polyps

Diverticula

Inflammatory bowel diseases (e.g., ulcerative colitis, Crohn disease)

Colonic stenosis secondary to ischemia, infection, or previous surgery

Perforated colon

Colonic fistula

Appendicitis

Extrinsic compression of the colon from extracolonic tumors (e.g., ovarian)

Extrinsic compression of the colon from an abscess

Malrotation of the gut

Colon volvulus

Intussusception

Hernia

### notes

## barium swallow

B

**Type of test** X-ray with contrast dye

### Normal findings

Normal size, contour, filling, patency, and positioning of the esophagus

### Test explanation and related physiology

The barium swallow provides a more thorough examination of the esophagus than most upper GI series (see p. 922). As in most barium contrast studies, defects in normal filling and narrowing of the barium column indicate tumor, strictures, or extrinsic compression from extraesophageal tumors or an abnormally enlarged heart and great vessels. Varices also can be seen as serpiginous, linear-filling defects. Such anatomic abnormalities as hiatal hernia, Schatzki rings, and diverticula (Zenker or epiphrenic) can be seen as well.

In patients with esophageal reflux, the radiologist may identify reflux of the barium from the stomach back into the esophagus. Muscular abnormalities (e.g., achalasia, diffuse esophageal spasm) can be detected easily by a barium swallow.

### Contraindications

- Patients with evidence of bowel obstruction
- Barium may create a stonelike impaction.
- *Patients with a perforated viscus or ruptured esophagus.*
  - If barium were to leak, the degree of infection would be much worse. Usually, when perforation is suspected, diatrizoate (Gastrografin), a water-soluble contrast medium, is used.
- Patients who are unable to cooperate for the test

### Potential complications

- Barium-induced fecal impaction

### Interfering factors

- Food within the esophagus prevents adequate visualization.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Instruct the patient not to take anything by mouth for at least 8 hours before testing. Usually the patient is kept NPO after midnight on the day of the test.

## 132 barium swallow

- Assess the patient's ability to swallow. If the patient tends to aspirate, inform the radiologist.

### During

- Note the following procedural steps:
  1. The fasting patient is asked to swallow the contrast medium. Usually this is barium sulfate in a milkshake-like substance.
  2. As the patient drinks the contrast through a straw, the x-ray table is tilted to the near-erect position.
  3. The patient is asked to roll into various positions so that the entire esophagus can be adequately visualized.
  4. With fluoroscopy, the radiologist follows the barium column through the entire esophagus.
- Note that this procedure is usually performed in the radiology department by a radiologist in approximately 30 minutes.

**PT** Tell the patient that no discomfort is associated with this test.

### After

**PT** Inform the patient of the need to evacuate all the barium. Cathartics are recommended. Initially stools are white but should return to a normal color with complete evacuation.

### Abnormal findings

Total or partial esophageal obstruction

Cancer

Scarred strictures

Lower esophageal rings

Peptic esophageal ulcers

Varices

Peptic or corrosive esophagitis

Achalasia

Esophageal motility disorders (e.g., presbyesophagus, diffuse esophageal spasm)

Diverticula

Chalasia

Extrinsic compression from extraesophageal tumors, cardiomegaly, or aortic aneurysm

### notes



**basic metabolic profile** (BMP, Chemistry panel, Chem 7, Sequential multiple analysis [SMA] 7, SMAC 7, Chemistry Screen)

B

The BMP is a frequently ordered panel of eight tests that provides information about the patient's metabolism, kidney function, blood glucose level, electrolytes, and acid-base balance. Depending on the reason for the test, blood may be drawn after fasting or on a random basis. Abnormal results can indicate the need for further testing. The tests are listed below and discussed separately:

**Electrolytes**

- Potassium (see p. 724)
- Chloride (see p. 233)
- Calcium (see p. 189)
- CO<sub>2</sub> (carbon dioxide, bicarbonate [HCO<sub>3</sub>], see p. 197)

**Kidney tests**

- BUN (see p. 155)
- Creatinine, serum (see p. 301)

**Glucose (see p. 462)**

Test results are often displayed as shown in the illustration below. For accuracy, check the normal ranges at the institution or lab where the test was performed.

<b>Na</b> 135-147 mEq/L	<b>Cl</b> 95-106 mEq/L	<b>BUN</b> 6-20 mg/dL	<b>Glucose</b> 70-110 mg/dL
<b>K</b> 3.5-5.2 mEq/L	<b>CO<sub>2</sub></b> 22-30 mEq/L	<b>Cr</b> 0.6-1.3 mg/dL	

**Bence Jones protein** (Free kappa and lambda light chains)

**Type of test** Urine

**Normal findings**

Kappa total light chain: < 0.68 mg/dL

Lambda total light chain: < 0.40 mg/dL

Kappa/lambda ratio: 0.7-6.2

**Test explanation and related physiology**

The detection of Bence Jones protein in the urine most commonly indicates multiple myeloma (especially when the urine levels are high). The test is used to detect and monitor the treatment and clinical course of multiple myeloma and other similar globulin diseases.

Bence Jones proteins are monoclonal light-chain portions of immunoglobulins found in 75% of patients with multiple myeloma. These proteins are made most notably by the plasma cells in these patients. They also may be associated with tumor metastases to the bone, chronic lymphocytic leukemias, lymphoma, macroglobulinemia, and amyloidosis.

Immunoglobulin light chains are usually cleared from the blood through the renal glomeruli and are reabsorbed in the proximal tubules; thus urine light-chain concentrations are normally very low or undetectable. The production of large amounts of monoclonal light chains, however, can overwhelm this reabsorption mechanism. Because the Bence Jones protein is rapidly cleared from the blood by the kidneys, it may be very difficult to detect in the blood; therefore urine is used for this study. Normally urine should contain no Bence Jones proteins.

Routine urine testing for proteins using reagent strips often does not reflect the type or amount of proteins in the urine. In fact, the strip may show a completely negative result despite large amounts of Bence Jones globulins in the urine. Proteins in the urine are best identified by *protein electrophoresis* of the urine and then are identified and quantified (*immunofixation*). Monitoring the urine M-spike (a spike on electrophoresis indicating multiple myeloma) is especially useful in patients with light-chain multiple myeloma in whom the serum M-spike may be very small or absent but in whom the urine M-spike is large.

**Interfering factors**

- Dilute urine may yield a false-negative result.
- ✚ High doses of aspirin or penicillin can cause false-positive results.

## Procedure and patient care

- See inside front cover for Routine Urine Testing.
- PT** Instruct the patient to collect an early morning specimen of at least 50 mL of uncontaminated urine in a container. It may be helpful to know the amount of these proteins excreted over 24 hours. If so, a 24-hour collection is ordered.
- Immediately transport the specimen to the laboratory. If it cannot be taken to the laboratory immediately, refrigerate it. Heat-coagulable proteins can decompose, causing a false-positive test.

## Abnormal findings

### ▲ Increased levels

- Multiple myeloma (plasmacytoma)
- Various metastatic tumors
- Chronic lymphocytic leukemia
- Amyloidosis
- Lymphoma
- Waldenström macroglobulinemia
- Osteogenic sarcoma
- Cryoglobulinemia
- Rheumatoid diseases

## notes

## 11 beta-prostaglandin F(2) alpha

**Type of test** Urine

### Normal findings

> 1000 ng/24 hours

### Test explanation and related physiology

Measurement of 11 beta-prostaglandin F(2) alpha in urine is useful in the evaluation of patients suspected of having systemic mastocytosis (systemic mast cell disease [SMCD]). SMCD is characterized by mast cell infiltration of extracutaneous organs (usually the bone marrow). Focal mast cell lesions in the bone marrow are found in approximately 90% of adult patients with SMCD.

Prostaglandin D(2) (PGD[2]) is generated by human mast cells, activated alveolar macrophages, and platelets. Although the most definitive test for SMCD is bone marrow biopsy (see p. 160), measurement of mast cell mediators such as beta prostaglandin in urine is advised for the initial evaluation of suspected cases. Elevated levels of 11 beta-prostaglandin F(2) alpha in urine are not specific for SMCD and may be found in patients with angioedema, diffuse urticaria, or myeloproliferative diseases in the absence of diffuse mast cell proliferation.

### Procedure and patient care

- See inside front cover for Routine Urine Testing.

### Abnormal findings

#### ▲ Increased levels

Systemic mast cell disease (SMCD)

notes

## bilirubin

**Type of test** Blood

### Normal findings

Adult/elderly/child:

Total bilirubin: 0.3-1.0 mg/dL or 5.1-17  $\mu\text{mol/L}$  (SI units)

Indirect bilirubin: 0.2-0.8 mg/dL or 3.4-12.0  $\mu\text{mol/L}$   
(SI units)

Direct bilirubin: 0.1-0.3 mg/dL or 1.7-5.1  $\mu\text{mol/L}$   
(SI units)

Newborn:

Total bilirubin: 1.0-12.0 mg/dL or 17.1-205  $\mu\text{mol/L}$   
(SI units)

### Possible critical values

Total bilirubin

Adult: > 12 mg/dL

Newborn: > 15 mg/dL


### Test explanation and related physiology

Bile, which is formed in the liver, has many constituents, including bile salts, phospholipids, cholesterol, bicarbonate, water, and bilirubin. Bilirubin metabolism begins with the breakdown of red blood cells (RBCs) in the reticuloendothelial system (Figure 5). Hemoglobin is released from RBCs and broken down to heme and globin molecules. Heme is then catabolized to form biliverdin, which is transformed into bilirubin. This form of bilirubin is called *unconjugated (indirect) bilirubin*. In the liver, indirect bilirubin is conjugated with a glucuronide, resulting in *conjugated (direct) bilirubin*. The conjugated bilirubin is then excreted from the liver cells and into the bile ducts and then into the bowel.

*Jaundice* is the discoloration of body tissues caused by abnormally high blood levels of bilirubin. This yellow discoloration is recognized when the total serum bilirubin exceeds 2.5 mg/dL.

*Physiologic jaundice of the newborn* occurs if the newborn's liver is immature and does not have enough conjugating enzymes.

When the jaundice is recognized either clinically or chemically, it is important (for therapy) to differentiate whether it is predominantly caused by unconjugated or conjugated bilirubin. This in turn will help differentiate the etiology of the defect. In general, jaundice caused by hepatocellular dysfunction (e.g., hepatitis) results in elevated levels of unconjugated bilirubin.

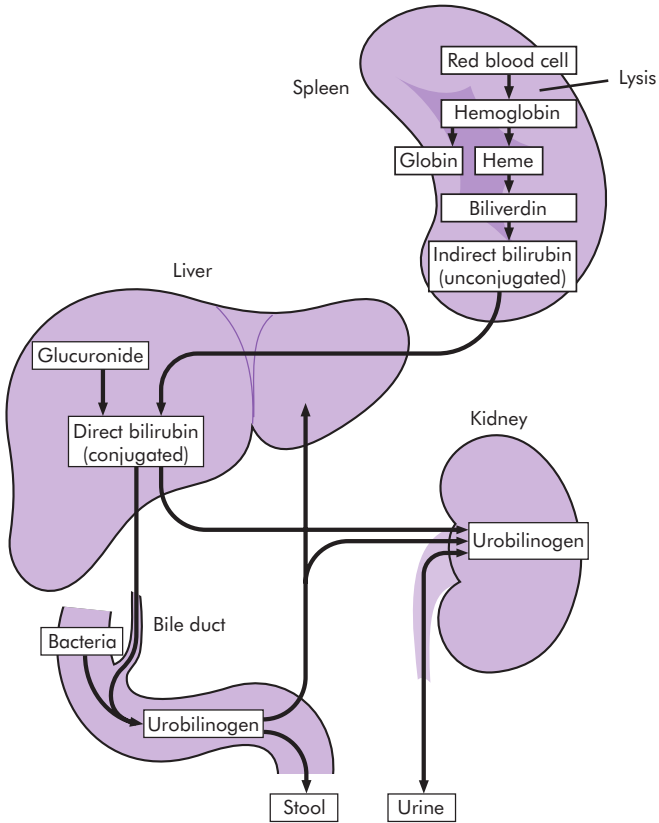

**FIG. 5** Bilirubin metabolism and excretion. The spleen, liver, kidneys, and gastrointestinal tract contribute to this process.

Jaundice resulting from extrahepatic obstruction of the bile ducts (e.g., gallstones or tumor blocking the bile ducts) results in elevated conjugated bilirubin levels; this type of jaundice usually can be resolved surgically or endoscopically.

The total serum bilirubin level is the sum of the conjugated (direct) and unconjugated (indirect) bilirubin. These are separated out when *fractionation* or *differentiation* of the total bilirubin to its direct and indirect parts is requested from the laboratory. Normally the unconjugated bilirubin makes up 70% to 85% of the total bilirubin. In patients with jaundice, when more

than 50% of the bilirubin is conjugated, it is considered a conjugated hyperbilirubinemia from gallstones, tumors, inflammation, scarring, or obstruction of the extrahepatic ducts. Unconjugated hyperbilirubinemia exists when less than 15% to 20% of the total bilirubin is conjugated. Diseases that typically cause this form of jaundice include accelerated erythrocyte (RBC) hemolysis or hepatitis.

### Interfering factors

- Blood hemolysis and lipemia can produce erroneous results.
- ✚ Drugs that may cause *increased* levels of total bilirubin include allopurinol, anabolic steroids, antibiotics, antimalarials, ascorbic acid, azathioprine, chlorpropamide, cholinergics, codeine, dextran, diuretics, epinephrine, meperidine, methotrexate, methyldopa, monoamine oxidase inhibitors, morphine, nicotinic acid (large doses), oral contraceptives, phenothiazines, quinidine, rifampin, salicylates, steroids, sulfonamides, theophylline, and vitamin A.
- ✚ Drugs that may cause *decreased* levels of total bilirubin include barbiturates, caffeine, penicillin, and salicylates (large doses).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: verify with laboratory
- Blood tube commonly used: red
- Use a heel puncture for blood collection in infants.
- Prevent hemolysis of blood during phlebotomy.
- Do *not* shake the tube; inaccurate test results may occur.
- Protect the blood sample from bright light. Prolonged exposure (> 1 hour) to sunlight or artificial light can reduce bilirubin content.

### Abnormal findings

#### ▲ Increased levels of conjugated (direct) bilirubin

Gallstones

Extrahepatic duct obstruction (tumor, inflammation, gallstone, scarring, or surgical trauma)

Extensive liver metastasis

Cholestasis from drugs

Dubin–Johnson syndrome

Rotor syndrome

## 140 bilirubin

### ▲ Increased levels of unconjugated (indirect) bilirubin

Hemolytic disease of the newborn

Hemolytic jaundice

Large-volume blood transfusion

Resolution of a large hematoma

Hepatitis

Sepsis

Neonatal hyperbilirubinemia

Hemolytic anemia

Crigler–Najjar syndrome

Gilbert syndrome

Pernicious anemia

Cirrhosis

Transfusion reaction

Sickle cell anemia

### notes



## bioterrorism infectious agents testing

B

**Type of test** Various (e.g., blood, urine, stool, tissue culture, sputum, lymph node biopsy, skin)

### Normal findings

Negative for evidence of infectious agent

### Test explanation and related physiology

There are many infectious agents used in bioterrorism, and it would be difficult to discuss each possible agent. In this test, those agents to which humans are most likely to be exposed, either in war or in a civilian terrorist attack, are discussed. Refer to [Table 4](#) for specific information on each agent. All documented cases must be reported to the Department of Public Health.

### Botulism infection

The botulinum toxin produced by *Clostridium botulinum* causes this disease. The gastrointestinal tract usually absorbs this organism after eating undercooked meat or sauces exposed to room temperature for prolonged periods. The organism also can be inhaled by handling these items or by open wound contamination of soil that contains *C. botulinum*.

The test used to diagnose this disease involves the identification of the toxin in the blood, stool, or vomitus of the affected individual. The food itself can also be tested.

### Anthrax

Anthrax is caused by *Bacillus anthracis*, which is a spore-forming gram-positive rod. Gastrointestinal anthrax is contracted by eating undercooked meat. Pulmonary anthrax results from inhalation of spores or tissues from infected animals. Once inhaled, it is always fatal without treatment. Cutaneous anthrax occurs after contact with contaminated meat, wool, hides, or leather from infected animals.

Culturing the organism in sheep blood agar makes the diagnosis. Appropriate specimens for culture would be stool, blood, sputum, and the cutaneous vesicle. Treatment for this disease is early institution of antibiotics and supportive care.

### Hemorrhagic fever (yellow fever)

This disease complex has many causative viruses, including arenavirus, bunyavirus (including hantavirus), filovirus (including Ebola), and flavivirus. The diagnosis is determined by clinical evaluation. However, viral cultures, serology, and immunohistochemistry of tissue specimens are possible. There is no specific treatment other than supportive medical therapy.

**TABLE 4** Bioterrorism infectious agents testing

Infection/infectious agent	Site of entry	Sources	Specimen	Tests
Botulism/ <i>Clostridium botulinum</i>	GI mucosal surfaces, lung, wound contamination	Undercooked meats, soil, dust	Blood, stool, vomitus, food	Botulinum toxin, mouse bioassay
Anthrax/ <i>Bacillus anthracis</i>	Lung, GI, skin	Undercooked meats, inhalation of spores from animal products, skin	Sputum, blood, stool, skin vesicle, food, spores	Culture, Gram stain
Yellow fever/Hantavirus, Ebola virus, multiple other viruses	Skin bite	Rodent or mosquito bites	Blood, sputum, tissue	Culture, serology for viral antigens
Plague infections/ <i>Yersinia pestis</i>	Skin bite	Infected fleas	Blood, sputum, lymph node aspirate	Culture of organism
Brucellosis/ <i>Brucella abortus</i> , <i>B. canis</i> , and so on	GI, lung, wound	Infected meat and milk products	Blood, sputum, food	Culture of organism
Smallpox/Variola virus	Lungs	Respiratory droplets, direct contact, contaminated clothing	Vesicle	Viral culture or viral identification with electron microscopy
Tularemia/ <i>Francisella tularensis</i>	Skin, GI tract, lungs	Ingestion of contaminated plants or water	Blood, sputum, stool	Culture of organism

## Plague

This disease is caused by *Yersinia pestis* and has three forms: bubonic (enlarged lymph nodes), septicemic (blood-borne), and pneumonic (aerosol). Pneumonic is by far the deadliest form of the infection.

The diagnosis is made by culture of the blood, sputum, or lymph node aspirate. This disease complex can be treated with antibiotics when started early in the course of the disease.

## Brucellosis

This disease is caused by *Brucella abortus*, *suis*, *melitensis*, or *canis*. It is contracted by ingestion of contaminated milk products, direct puncture of the skin (in butchers and farmers), or by inhalation. *Brucella* can be cultured from a blood, sputum, or food specimen. Serology testing is also possible.

## Smallpox

Smallpox is a serious, contagious, and sometimes fatal infectious disease caused by the variola virus (a DNA virus). There is no specific treatment for smallpox disease, and the only prevention is vaccination. It is very easily spread and is therefore considered a potential bioterrorism weapon.

Viral culture, serology, immunohistochemistry, or electron microscopy can make the diagnosis. The best specimen is the vesicular rash. There is no treatment for the disease, but vaccination is available and is offered to all those at risk for bioterrorism.

## Tularemia

This disease is caused by a bacterium called *Francisella tularensis*. It is contracted by drinking contaminated water or eating vegetation contaminated by infected animals. Diagnosis is made by culture of the blood, sputum, or stool.

## Procedure and patient care

### Before

- Maintain strict adherence to all procedures to avoid violations in isolation or contamination.
- Laboratory personnel must strictly adhere to all universal standard precautions and biohazard concerns.

### During

- If an enema is used to obtain a botulinum stool specimen, use sterile water. Saline can negate results.
- Send enough blood for adequate testing. Usually two red-top tubes are adequate. It is best to send it on ice.
- Send food for testing in its original container.

## 144 bioterrorism infectious agents testing

- For anthrax or smallpox testing of a cutaneous lesion, soak one or two culture swabs with fluid from an unopened lesion.

### After

- Identify all potential sources of contamination.
- Isolate individuals suspected of having a contagious disease.

### Abnormal findings

See [Table 4](#).

**bladder cancer markers** (Bladder tumor antigen [BTA], Nuclear matrix protein 22 [NMP22])

**Type of test** Urine

**Normal findings**

BTA: < 14 units/mL

NMP22: < 10 units/mL

FISH: No chromosomal amplification or deletions noted

**Test explanation and related physiology**

The recurrence rate for superficial bladder cancers that have been resected by transurethral cystoscopy is high. Surveillance testing requires frequent urine testing for cytology and frequent cystoscopic evaluations. The use of bladder tumor markers may provide an easier, cheaper, and more accurate method of diagnosing recurrent bladder cancer.

Bladder tumor antigen (BTA) and nuclear matrix protein 22 (NMP22) are proteins produced by bladder tumor cells and deposited into the urine. Normally none or very low levels of these proteins are found in the urine. When levels of bladder cancer tumor markers are normal, cystoscopy rarely yields positive results. When these markers are elevated, bladder tumor recurrence is strongly suspected and cystoscopy is indicated to confirm bladder cancer recurrence.

NMP22 may also be a good screening test for patients at increased risk for developing bladder cancer. However, these markers can be elevated in other circumstances (i.e., recent urologic surgery, urinary tract infection, calculi). Cancers involving the ureters and renal pelvis may also be associated with increased BTA and NMP22.

Bladder cancer cells have been found to exhibit aneuploidy (gene amplifications on chromosomes 3, 7, and 17, and the loss of the 9p21 locus on chromosome 9). Using DNA probes, through *fluorescence in situ hybridization* (FISH), these chromosomal abnormalities can be identified with great accuracy. FISH can be performed on cells isolated in a fresh urine specimen or cells available on a ThinPrep slide (similar to Pap smears [see p. 667]).

Although not a tumor marker, a cytology test is available that can be used in the early detection of bladder cancer recurrence. It is an immunocytofluorescence technique based on monoclonal antibodies to two antigens: a mucin glycoprotein and a carcinoembryonic antigen (CEA). These antigens are expressed by

## 146 bladder cancer markers

tumor cells found in patients with bladder cancer and are exfoliated in the urine.

### Interfering factors

- These proteins are very unstable. If the urine is not immediately stabilized, false negatives may occur.
- Active infection (including sexually transmitted diseases) of the lower urologic tract can cause false elevations.
- Kidney or bladder calculi can cause false elevations.

### Procedure and patient care

- See inside front cover for Routine Urine Testing.

**PT** Tell the patient that no fasting is required.

- A urine specimen should be collected, preferably from the first void of the day.
- The specimen should be transported to the laboratory immediately to avoid deterioration of the cells.
- If a time delay is required, the specimen should be refrigerated.

### Abnormal findings

Bladder cancer

Nonbladder urologic cancer (e.g., ureters, renal pelvis)

notes

## blood culture and sensitivity

B

**Type of test** Blood

### Normal findings

Negative

### Test explanation and related physiology

Blood cultures are obtained to detect the presence of bacteria in the blood (bacteremia). An episode of bacteremia is usually accompanied by chills and fever; thus the blood culture should be drawn when the patient manifests these signs to increase the chances of growing bacteria on the cultures. It is important that at least two culture specimens be obtained from two different sites. If one produces bacteria and the other does not, it is safe to assume that the bacteria in the first culture may be a contaminant and not the infecting agent. When both cultures grow the infecting agent, bacteremia exists and is caused by the organism that is growing in the culture.

If the patient is receiving antibiotics during the time that the cultures are drawn, the laboratory should be notified. Resin may be added to the culture medium to negate the antibiotic effect in inhibiting growth of the offending bacteria in the culture. If cultures are to be performed while the patient is on antibiotics, the blood culture specimen should be taken shortly before the next dose of the antibiotic is administered. All cultures preferably should be performed before antibiotic therapy is initiated.

Culture specimens drawn through an IV catheter are frequently contaminated, and tests using them should not be performed unless catheter sepsis is suspected. Most organisms require approximately 24 hours to grow in the laboratory, and a preliminary report can be given at that time. Often, 48 to 72 hours is required for growth and identification of the organism. Anaerobic organisms may take longer to grow.

### Interfering factors

- Contamination of the blood specimen, especially by skin bacteria, may occur.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

**PT** Tell the patient that no fasting is required.

## 148 blood culture and sensitivity

### During

- Carefully prepare the venipuncture site with povidone-iodine (Betadine). Allow the skin to dry.
- Clean the tops of the Vacutainer tubes or culture bottles with povidone-iodine and allow them to dry. Some laboratories suggest cleaning with 70% alcohol after cleaning with povidone-iodine and air drying.
- Collect approximately 10 to 15 mL of venous blood by venipuncture from each site in a 20-mL syringe.
- Discard the needle on the syringe and replace it with a second sterile needle before injecting the blood sample into the culture bottle.
- Inoculate the anaerobic bottle first if both anaerobic and aerobic cultures are needed.
- Mix gently after inoculation.
- Label the specimen with the patient's name, date, time, and tentative diagnosis.
- Indicate on the laboratory slip any medications that may affect test results.

### After

- Transport the culture bottles immediately to the laboratory (or at least within 30 minutes).
- Notify the physician of any positive results so that appropriate antibiotic therapy can be initiated.

### Abnormal findings

Bacteremia

notes



**blood smear** (Peripheral blood smear, Red blood cell [RBC] morphology, RBC smear)

**Type of test** Blood

**Normal findings**

Normal quantity of red and white blood cells (RBCs, WBCs) and platelets

Normal size, shape, and color of RBCs

Normal WBC differential count

Normal size and granulation of platelets

**Test explanation and related physiology**

Examination of the peripheral blood smear can provide a significant amount of information concerning drugs and diseases that affect erythrocytes (RBCs), leukocytes (WBCs), or platelets. Furthermore, other congenital and acquired diseases can be diagnosed. When special stains are applied to the blood smear, leukemia, infection, infestation, and other diseases can be identified.

When adequately prepared and examined microscopically by an experienced technologist or pathologist, a smear of peripheral blood can be very diagnostic. All three hematologic cell lines—RBCs, platelets, and WBCs—can be examined. In the peripheral blood, five different types of WBCs can routinely be identified: neutrophils, eosinophils, basophils, lymphocytes, and monocytes. The first three are also referred to as *granulocytes*. Please see the discussion in bone marrow biopsy (see p. 160) for more information concerning the various elements of blood.

Microscopic examination of the RBCs can reveal variations in RBC size (anisocytosis), shape (poikilocytosis), color, or intracellular content (Box 3). Classification of RBCs according to these variables is most helpful in identifying the causes of anemia and the presence of other diseases.

The WBCs are examined for total quantity, differential count, and degree of maturity. An increased number of immature WBCs may indicate leukemia or infection. A decreased WBC count indicates a failure of marrow to produce WBCs (because of drugs, chronic disease, neoplasia, or fibrosis), peripheral destruction, or sequestration.

**Platelet examination**

Finally, an experienced laboratory technologist also can estimate platelet number. Platelets are small cell fragments that do not contain a nucleus. The contents of the granules in a platelet are released to promote clotting.

**BOX 3 Microscopic examination of red blood cells**

**RBC size abnormalities**

Microcytes (small RBCs)

- Iron deficiency
- Thalassemia
- Hemoglobinopathies

Macrocytes (larger size)

- Vitamin B<sub>12</sub> or folic acid deficiency
- Reticulocytosis secondary to increased erythropoiesis (RBC production)
- Occasional liver disorder

**RBC shape abnormalities**

Spherocytes (small and round)

- Hereditary spherocytosis
- Acquired immunohemolytic anemia

Elliptocytes (crescent)

- Iron deficiency
- Hereditary elliptocytosis

Codocytes or target cells (thin cells with less hemoglobin)

- Hemoglobinopathies
- Thalassemia

Echinocytes (Burr cells)

- Uremia
- Liver disease

**RBC color abnormalities**

Hypochromic (pale)

- Iron deficiency
- Thalassemia

Hyperchromasia (more colored)

- Concentrated hemoglobin, usually caused by dehydration

**RBC intracellular structure**

Nucleated (normoblasts). (Mature RBCs are round with a small central pallor without any intracellular structures. They do not have a nucleus. Immature RBCs [reticulocytes] do contain intracellular RNA. Immature nucleated cells are not normally found in the peripheral blood and indicate increased RBC synthesis.)

- Anemia
- Chronic hypoxemia
- Normal for an infant
- Marrow-occupying neoplasm or fibrotic tissue

Basophilic stippling (refers to bodies enclosed or included in the cytoplasm of the RBCs)

- Lead poisoning
- Reticulocytosis

Howell–Jolly bodies (small, round remnants of nuclear material remaining within the RBC)

- After a surgical splenectomy
- Hemolytic anemia
- Megaloblastic anemia
- Functional asplenia (after splenic infarction)

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting is required.

### During

- Collect a drop of blood from a finger stick or heel stick and place it on a slide.
- If necessary, perform a venipuncture and collect the blood in a lavender-top tube.
- Note that a blood smear is first studied with an automated cytometer programmed to recognize abnormal blood cell shapes and other variations. An evaluation smear is performed by a technologist. Low counts may be hand counted to ensure accuracy. The most accurate smear requires review by a pathologist.

### After

- Apply pressure to the venipuncture site.

### Abnormal findings

See [Box 3](#) in the Test explanation and related physiology section.

### notes

**blood typing** (Blood group microarray testing)**Type of test** Blood**Normal findings**

Compatibility

**Test explanation and related physiology**

With blood typing, ABO and Rh antigens can be detected in the blood of prospective blood donors and potential blood recipients. This test is also used to determine the blood type of expectant mothers and newborns. A description of the ABO system, Rh factors, and blood crossmatching is reviewed here.

**ABO system**

Human blood is grouped according to the presence or absence of A or B antigens. The surface membranes of group A red blood cells (RBCs) contain A antigens; group B RBCs contain B antigens on their surface; group AB RBCs have both A and B antigens; and group O RBCs have neither A nor B antigens. In general, a person's serum does not contain antibodies to match the surface antigen on their RBCs. That is, persons with group A antigens (type A blood) will not have anti-A antibodies; however, they will have anti-B antibodies. The converse is true for persons with group B antigens. Group O blood will have both anti-A and anti-B antibodies (Table 5). These antibodies against A and B blood group antigens are formed in the first 3 months of life after exposure to similar antigens on the surface of naturally occurring bacteria in the intestine.

**TABLE 5** Blood typing

Blood type (ABO, Rh)	Antigens present	Antibodies possibly present	Percent of general population
O, +	Rh	A, B	35
O, - <sup>a</sup>	None	A, B, Rh	7
A, +	A, Rh	B	35
A, -	A	B, Rh	7
B, +	B, Rh	A	8
B, -	B	A, Rh	2
AB, + <sup>b</sup>	A, B, Rh	None	4
AB, -	A, B	Rh	2

<sup>a</sup>Universal donor.<sup>b</sup>Universal recipient.

Blood transfusions are actually transplantations of tissue (blood) from one person to another. It is important that the recipient not have antibodies to the donor's RBCs. If this were to occur, there could be a hypersensitivity reaction, which can vary from mild fever to anaphylaxis with severe intravascular hemolysis. If donor ABO antibodies are present against the recipient antigens, usually only minimal reactions occur.

Persons with group O blood are considered *universal donors* because they do not have antigens on their RBCs. People with group AB blood are considered *universal recipients* because they have no antibodies to react to the transfused blood. Group O blood is usually transfused in emergent situations in which rapid, life-threatening blood loss occurs and immediate transfusion is required. The chance of a transfusion reaction is least when type O is used. Women of childbearing potential should receive group O-negative blood, and men generally receive group O-positive blood when emergency transfusion before type-specific or cross-matched blood is required.

ABO typing is not required for autotransfusions (blood donated by a patient several weeks before a major operation and then transfused postoperatively). However, in most hospitals, ABO typing is performed on those patients in the event that further blood transfusion of banked blood is required.

### Rh factor

The presence or absence of Rh antigens on the RBC's surface determines the classification of Rh positive or Rh negative. After ABO compatibility, Rh factor is the next most important antigen affecting the success of a blood transfusion. The major Rh factor is Rh<sub>0</sub>(D). There are several minor Rh factors. If Rh<sub>0</sub>(D) is absent, the minor Rh antigens are tested. If negative, the patient is considered *Rh negative (Rh-)*.

Rh- persons may develop antibodies to Rh antigens if exposed to Rh-positive (Rh+) blood by transfusions or fetal-maternal blood mixing. All women who are pregnant should have a blood typing and Rh factor determination. If the mother's blood is Rh-, the father's blood should also be typed. If his blood is Rh+, the woman's blood should be examined for the presence of Rh antibodies (by the indirect Coombs test; see p. 287). Hemolytic disease of the newborn can be prevented by Rh typing during pregnancy. If the mother is Rh-, she should be advised that she is a candidate for RhoGAM (Rh immunoglobulin that "neutralizes" the Rh antigen) after the birth. RhoGAM can reduce the chance of fetal hemolytic problems during subsequent pregnancies.

## 154 blood typing

### Other blood typing systems

There are nine different gene codes for blood groups assayed. Most are minor and not clinically significant. However, in certain clinical circumstances, these minor blood group antigens and acquired antigens can become significant. This may occur with frequent blood transfusions or in patients with leukemia or lymphoma.

### Blood crossmatching

Although typing for the major ABO and Rh antigens is no guarantee that a reaction will not occur, it does greatly reduce the possibility of such a reaction. Many potential minor antigens are not routinely detected during blood typing. If allowed to go unrecognized, these minor antigens also can initiate a blood transfusion reaction. Therefore blood is not only typed but also crossmatched to identify a mismatch of blood caused by minor antigens. Crossmatching always includes an indirect Coombs test. Only blood products containing RBCs need to be crossmatched. Plasma products do not need to be crossmatched but should be ABO compatible because other cells (WBCs and platelets) have ABO antigens.

Homologous (donor and recipient are different people) and directed (recipient chooses the donor) blood for donation must be rigorously tested before transfusion. Autologous (recipient and donor is the same person) blood for transfusions, however, is not subject to that same testing.

Finally, one must be aware of graft-versus-host disease (GVHD) in which donor lymphocytes included in the blood transfusion may engraft and multiply in the recipient.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red. Verify with laboratory.

### Abnormal findings

See [Test explanation and related physiology](#) section (see p. 153).

### notes

## blood urea nitrogen (BUN)

**Type of test** Blood

### Normal findings

Adult: 10-20 mg/dL or 3.6-7.1 mmol/L (SI units)

Elderly: may be slightly higher than those of adult

Child: 5-18 mg/dL

Infant: 5-18 mg/dL

Newborn: 3-12 mg/dL

Cord: 21-40 mg/dL

### Possible critical values

> 100 mg/dL (indicates serious impairment of renal function)

### Test explanation and related physiology

The BUN measures the amount of urea nitrogen in the blood. Urea is formed in the liver as the end product of protein metabolism. During ingestion, protein is broken down into amino acids. In the liver, these amino acids are catabolized, and free ammonia is formed. The ammonia is combined to form urea, which is then deposited into the blood and transported to the kidneys for excretion. Therefore BUN is directly related to the metabolic function of the liver and the excretory function of the kidney. It serves as an index of the function of these organs. Patients who have elevated BUN levels are said to have azotemia.

Nearly all renal diseases cause inadequate excretion of urea, which causes the blood concentration to rise above normal. BUN also increases in conditions other than primary renal disease. For example, when excess amounts of protein are available for hepatic catabolism (from gastrointestinal [GI] bleeding), large quantities of urea are made.

BUN is interpreted in conjunction with the creatinine test (see p. 301). These tests are referred to as *renal function studies*.

The BUN/creatinine ratio is a good measurement of kidney and liver function. The normal adult range is 6 to 25, with 15.5 being the optimal adult value for this ratio.

### Interfering factors

- Changes in protein intake may affect BUN levels.
- Advanced pregnancy may cause increased BUN levels.
- Overhydration and underhydration will affect BUN levels.
- GI bleeding can cause increased BUN levels.

## 156 blood urea nitrogen (BUN)

- ✦ Drugs that may cause *increased* BUN levels include allopurinol, aminoglycosides, cephalosporins, chloral hydrate, cisplatin, furosemide, guanethidine, indomethacin, methotrexate, methyl dopa, nephrotoxic drugs (e.g., amphotericin B, aspirin, bacitracin, carbamazepine, colistin, gentamicin, methicillin, neomycin, penicillamine, polymyxin B, probenecid, vancomycin), propranolol, rifampin, spironolactone, tetracyclines, thiazide diuretics, and triamterene.
- ✦ Drugs that may cause *decreased* BUN levels include chloramphenicol and streptomycin.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

##### Prerenal causes

Hypovolemia  
Shock  
Burns  
Dehydration  
Congestive heart failure  
Myocardial infarction  
GI bleeding  
Excessive protein ingestion  
Alimentary tube feeding  
Excessive protein catabolism  
Starvation  
Sepsis

##### Renal causes

Renal disease (e.g., glomerulonephritis, pyelonephritis, acute tubular necrosis)  
Renal failure  
Nephrotoxic drugs

##### Postrenal azotemia

Ureteral obstruction  
Bladder outlet obstruction

#### ▼ Decreased levels

Liver failure  
Overhydration caused by fluid overload or syndrome of inappropriate antidiuretic hormone (SIADH)  
Negative nitrogen balance (e.g., malnutrition or malabsorption)  
Pregnancy  
Nephrotic syndrome



## bone densitometry (Bone mineral content [BMC], Bone mineral density [BMD], DEXA scan)

**Type of test** X-ray

### Normal findings

Normal: < 1 standard deviation below normal ( $> -1$ )

Osteopenia: 1 to 2.5 standard deviations below normal ( $-1$  to  $-2.5$ )

Osteoporosis: > 2.5 standard deviations below normal ( $< -2.5$ )

### Test explanation and related physiology

Bone densitometry is used to determine bone mineral content (BMC) and density (BMD) to diagnose osteoporosis as early as possible. It is also used to monitor patients who are undergoing treatment for osteoporosis. *Osteoporosis* and *low bone mass* are terms used for bones that become weakened and fracture easily. This most commonly occurs in postmenopausal women. However, other diseases are associated with osteoporosis, such as malnourishment and osteopenic endocrinopathies (e.g., hyperparathyroidism). Bone densitometry can provide early and accurate measurements of bone strength based on BMD.

Several groups of bones are routinely evaluated. The lumbar spine is one of the best representatives of cancellous bone. The radius is the most frequently studied cortical bone. The proximal hip (neck of the femur) is the best representative of mixed (cancellous and cortical) bone. However, specific bone sites can be evaluated if they are particularly symptomatic.

*Dual-energy x-ray absorptiometry (DEXA)* is the method most commonly used. Because DEXA uses two photons, more energy is produced so that bones (spine and hip [femoral neck]) surrounded by a lot of soft tissue can be more easily penetrated. The source of the photon is placed on one side of the bone to be studied. The gamma detector is placed on the other side. Increased bone density is associated with increased bone photon absorption and therefore less photon recognition at the site of the gamma detector.

Several other methods are available to measure BMD. *Quantitative computed tomography (QCT)* uses CT technology to measure central bones, especially the spine. *Ultrasound absorption (quantitative ultrasound)* can be used to measure peripheral bones (heel [calcaneus], patella, or midtibia).

BMD is usually reported in terms of standard deviation (SD) from mean values. T scores compare the patient's results to a group of young healthy adults. Z scores compare the patient's

## 158 bone densitometry

results to a group of age-matched controls. The World Health Organization has defined *low bone mass* as a BMD value greater than 1 SD below peak bone mass levels in young women and *osteoporosis* as a value greater than 2.5 SD below that same measurement scale. Positive T scores indicate a normal BMD. Negative T scores indicate reduced BMD.

*Vertebral fracture assessment (VFA)* can be performed utilizing the images generated by the DEXA scan. Images of the lower thoracic and lumbar spine are examined. If a vertebral fracture is identified, bone mineral strengthening medications are recommended despite T score. Presence of a vertebral fracture indicates a substantial risk for a subsequent vertebral or nonvertebral fracture independent of the bone mineral density or other osteoporosis risk factors.

BMD testing is an important part of routine screening testing for postmenopausal women. In general, BMD is recommended every 2 years to screen for osteoporosis. Women and men with known osteoporotic fractures, hyperparathyroidism, or administration of long-term steroid therapy may benefit from annual BMD testing.

### Interfering factors

- Barium may falsely increase the density of the lumbar spine. BMD measurements should not be performed for about 10 days after barium studies.
- Calcified abdominal aortic aneurysm may falsely increase BMD of the spine.
- Internal fixation devices of the hip or radius will falsely increase BMD of those bones.
- Overlying metal jewelry or other objects may falsely increase BMD.
- Previous fractures or severe arthritis changes of the bone can falsely increase BMD.
- Metallic clips placed in the vertebrae of patients who have had previous abdominal surgery can falsely increase BMD.
- Prior bone scans can falsely decrease BMD because the photons generated from the bone (as a result of the previously administered bone scan radionuclide) will be detected by the scintillator detector.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.

- PT** Tell the patient that no fasting or sedation is required.
- PT** Instruct the patient to remove all metallic objects (e.g., belt buckles, zippers, coins, keys) that might be in the scanning path.

### During

- Note the following procedural steps:
  1. The patient lies supine on an imaging table, with his or her legs supported and placed on a padded box to flatten the pelvis and lumbar spine.
  2. Under the table, a photon generator is slowly and successively passed under the lumbar spine.
  3. A scintillator (gamma or x-ray) detector/camera is passed over the patient in a manner parallel to that of the generator. An image of the lumbar spine and hip bone is obtained by the scintillator camera and projected onto a computer monitor.
  4. Next, the appropriate foot is applied to a brace that internally rotates the nondominant hip, and the procedure is repeated over the hip. A similar procedure is performed for radius evaluation.
  5. When the radius is examined, the nondominant arm is preferred unless there is a history of fracture to that bone.
- Note that the data are interpreted and reported by a radiologist or a physician trained in nuclear medicine.
- Note that BMD studies take about 30 minutes to perform and are free of any discomfort. Only minimal radiation is used for this procedure.
- Note that there are numerous types of bone densitometry machines. Peripheral units that quickly scan the finger, heel, or forearm are often used to identify patients at risk for osteoporosis. Abnormal results are followed up with the more comprehensive table procedure previously described.

### After

- On the computer screen, a small window of the lumbar spine, femoral neck, or distal radius is drawn. The computer calculates the amount of photons not absorbed by the bone. This is the BMC. BMD is computed as follows:

$$\text{BMD} = \frac{\text{BMC}(\text{g}/\text{cm}^2)}{\text{Surface area of the bone}}$$

### Abnormal findings

Low bone mass  
Osteoporosis

### notes

## bone marrow biopsy (Bone marrow examination, Bone marrow aspiration)

**Type of test** Microscopic examination of tissue

### Normal findings

Active erythroid cell line, myeloid and lymphoid cell lines, and megakaryocyte (platelet) production. See range of cell types below.

Cell type	Range (%)
Neutrophilic series	49.2-65.0
Myeloblasts	0.2-1.5
Promyelocytes	2.1-4.1
Myelocytes	8.2-15.7
Eosinophilic series	1.2-5.3
Myelocytes	0.2-1.3
Metamyelocytes	0.4-2.2
Bands	0.2-2.4
Segmented	0.0-1.3
Basophilic and mast cells	0.0-0.2
Erythrocyte series	18.4-33.8
Pronormoblasts	0.2-1.3
Basophilic	0.5-2.4
Polychromatophilic	17.9-29.2
Orthochromatic	0.4-4.6
Monocytes	0.0-0.8
Lymphocytes	11.1-23.2
Plasma cells	0.4-3.9
Megakaryocytes	0.0-0.4
Reticulum cells	0.0-0.9
Myeloid to erythrocyte (M/E) ratio	2:1-4:1

Normal iron content is demonstrated by staining with Prussian blue.

### Test explanation and related physiology

Bone marrow examination is an important part of the evaluation of patients with hematologic diseases. Indications for bone marrow examination include the following:

- To evaluate anemias, leukopenia, or thrombocytopenia
- To diagnose leukemia, myelodysplastic syndromes, myeloproliferative neoplasms, and plasma cell dyscrasia

- To document abnormal iron stores
- To document bone marrow infiltrative diseases (e.g., neoplasm, infection, or fibrosis)
- To stage lymphomas or other cancers

The bone marrow is located in the central fatty core of cancellous bone (particularly sternum, rib, and pelvis). There, the blood-forming cells produce blood cells and release them into the circulation.

Examination of bone marrow reveals the number, size, and shape of the RBCs, WBCs, and megakaryocytes (platelet precursors) as these cells evolve through their various stages of development in the bone marrow. Samples of bone marrow can be obtained by aspiration, bone marrow biopsy, or surgical removal. An aspiration is obtained for cell morphology, immunophenotyping, cytogenetics, or microbiology cultures. Microscopic examination of the marrow biopsy includes estimation of cellularity, identification of disordered hematopoiesis, and determination of the presence of infiltrative diseases (fibrosis or neoplasms, both primary and metastatic). Estimation of iron storage is performed on bone marrow aspirates or nondecalcified clot sections.

For the estimation of cellularity, the specimen is examined and the relative quantity of each cell type is determined. Leukemias or leukemoid drug reactions are suspected when increased numbers of leukocyte precursors are present. Physiologic marrow leukemoid compensation is also seen with infection. Decreased numbers of marrow leukocyte precursors occur in patients with myelofibrosis, metastatic neoplasia, or agranulocytosis/aplastic anemia; in elderly patients; and after radiation therapy or chemotherapy. Some drugs or infections can diminish leukocyte production.

Increased numbers of marrow RBC precursors occur with polycythemia vera or as physiologic compensation to blood loss (hemorrhage or hemolysis). Decreased numbers of marrow RBC precursors occur with erythroid hypoplasia after chemotherapy, infection (parvovirus), aplastic anemia, radiation therapy, administration of other toxic drugs, iron administration, or marrow replacement by fibrotic tissue or neoplasms.

Increased numbers of platelet precursors (megakaryocytes) can be the result of compensation to platelet loss from a recent hemorrhage. They are also seen in some forms of acute and chronic myeloid leukemias. This increase also may be compensatory in patients with platelet sequestration (secondary hypersplenism associated with portal hypertension) or platelet destruction (idiopathic thrombocytopenic purpura). Platelet

## 162 bone marrow biopsy

counts decrease, and the marrow compensates by increasing production. Decreased numbers of megakaryocytes occur in patients who have had radiation therapy, chemotherapy, or other drug therapy and in patients with neoplastic or fibrotic marrow infiltrative diseases. Patients with aplastic anemia also have decreased numbers of megakaryocytes.

Increased numbers of lymphocyte precursors occur in chronic, viral, or mycoplasmal infections (e.g., mononucleosis), lymphocytic leukemia, and lymphoma. Plasma cells and lymphocytes are increased in patients with plasma cell dyscrasia, lymphomas, hypersensitivity states, autoimmune disease, chronic infections, and other chronic inflammatory diseases.

Estimation of cellularity also can be expressed as a ratio of myeloid precursors to erythroid precursors (M/E ratio). The normal M/E ratio is approximately 3:1. The M/E ratio is greater than normal in those diseases in which leukocyte precursors are increased or erythroid precursors are decreased. The M/E ratio is below normal when either leukocyte precursors are decreased or erythroid precursors are increased.

Bone marrow aspiration and biopsy are performed by a physician or midlevel health-care provider. The duration of these procedures is approximately 20 minutes. The patient may have some apprehension when pressure is applied to puncture the outer table of the bone during biopsy-specimen removal or aspiration. The patient probably will feel pain during lidocaine infiltration and pressure when the syringe plunger is withdrawn for aspiration.

### Contraindications

- Patients with acute coagulation disorders, because of the risk of excessive bleeding
- Patients who cannot cooperate or remain still during the procedure

### Potential complications

- Hemorrhage, especially if the patient has a coagulopathy
- Infection, especially if the patient is leukopenic

### Procedure and patient care

#### Before


**PT** Explain the procedure to the patient.

- Obtain a written informed consent for this procedure.
- Encourage the patient to verbalize fears because many patients are anxious concerning this study. Conscious sedation may be required.

- Assess the results of the coagulation studies. Report any evidence of coagulopathy to the physician. Platelets should be greater than 20,000, and international normalized ratio (INR) should be less than 1.5.
  - Obtain an order for sedatives if the patient appears extremely apprehensive.
- PT** Instruct the patient to remain very still throughout the procedure.

### During

- PT** Inform the patient that during bone marrow aspiration, most patients feel pain or a burning sensation during lidocaine infiltration and pressure when the syringe plunger is withdrawn for aspiration.
- Conscious sedation may be provided for this procedure.
  - Note the following procedural steps for *bone marrow aspiration*:
    1. The procedure is usually begun as described in step 1 below for bone marrow biopsy.
    2. For aspiration, an Illinois- or Jamshidi-type large-bore needle is used.
    3. When inside the marrow, a syringe is used to aspirate marrow contents (Figure 6).
    4. Several small-volume (0.5-2 mL) samples of bone marrow are aspirated.
    5. The aspirate is placed in an appropriate blood specimen collecting test tube depending on the test requested.
  - Note the following procedural steps for *bone marrow biopsy*:

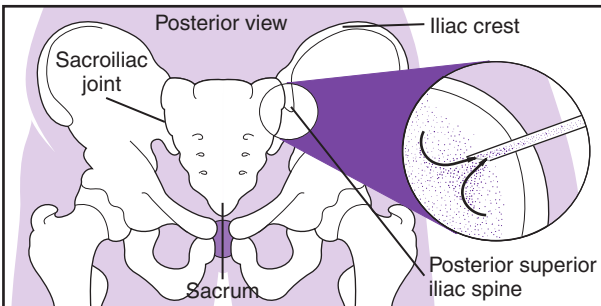

**FIG. 6** Bone marrow aspiration. Samples of the bone marrow are taken from along the posterior superior iliac spine.

## 164 bone marrow biopsy

1. The skin and soft tissues overlying the posterior superior spine of the iliac bone are prepped and draped. A small skin incision is made in that area after local anesthesia is provided.
2. A Jamshidi (or other specialized needle) needle is positioned into the bone.
3. The aspiration specimen is obtained first. With repositioning of the needle (to avoid aspiration artifact), the biopsy specimen is obtained and placed in a formalin fixative. It is then sent to the pathology laboratory for analysis.
4. Bilateral bone marrow biopsies may be performed for staging of lymphoma or other neoplasms.

### After

- Apply pressure to the puncture site to arrest any bleeding.
- Apply an adhesive bandage.
- Observe the puncture site for bleeding. Ice packs may be used to minimize bleeding.
- Assess for tenderness and erythema, which may indicate infection. Report this to the physician.
- Normally, place the patient in the supine position at bed rest for 30 to 60 minutes after the test. This provides pressure on the biopsy site.
- Note that some patients complain of tenderness at the puncture site for several days after this study. Mild analgesics may be ordered.

### Abnormal findings

Metastatic neoplasm	Lymphoblastic leukemia
Myeloproliferative neoplasms	Myeloid leukemia
Infection (e.g., viral, bacterial, fungal)	Myelodysplastic syndromes
Agranulocytosis	Anemia
Red cell aplasia	Chronic inflammatory disease
Plasma cell dyscrasia	Acquired immunodeficiency syndrome (AIDS)
Hodgkin lymphoma	Bone marrow aplasia
Non-Hodgkin lymphoma	

### notes



## bone scan

**Type of test** Nuclear scan

### Normal findings

No evidence of abnormality

### Test explanation and related physiology

The bone scan permits examination of the skeleton by a scanning camera after IV injection of a radionuclide material. After injection, the radiopharmaceutical is taken up by the bone. A realistic image of the bones is apparent.

The degree of radionuclide uptake is related to the metabolism of the bone. Normally a uniform concentration should be seen throughout the bones of the body. There is symmetric distribution of activity throughout the skeletal system in healthy adults. An increased uptake of radionuclide is abnormal and may represent tumor, arthritis, fracture, degenerative bone and joint changes, osteomyelitis, bone necrosis, osteodystrophy, and Paget disease. These areas of concentrated radionuclide uptake are often called *hot spots* and are detectable months before an ordinary x-ray image can reveal the pathology. Hot spots occur because new bone growth is usually stimulated around areas of pathology. If pathology exists and there is no new bone formation around the lesion, the scan will not pick up the abnormality. Increased uptake of radionuclide is also seen in the normal physiologic active epiphyses of children.

The major reason a bone scan is performed is to detect metastatic cancer to the bone. Bone scans are also useful in staging primary bone tumors. Bone scans may be serially repeated to monitor tumor response to antineoplastic therapy.

Bone scans also provide valuable information in the evaluation of patients with trauma or unexplained pain. Bone scanning is much more sensitive than routine x-ray images in detecting small and difficult-to-find fractures, especially in the spine, ribs, face, and small bones of the extremities. Bone scans are used to determine the age of a fracture as well.

Although the bone scan is extremely sensitive, unfortunately it is not very specific. Fractures, infections, tumors, and arthritic changes all appear similar in this scan. A three-phase bone scan may be performed if inflammation (arthritis) or infection (osteomyelitis, septic arthritis) is suspected. In a three-phase bone scan, imaging is performed at three different times after injection of the radionuclide. Early uptake of the radionuclide would indicate

infection or inflammation rather than neoplasm. Uptake of the radionuclide on delayed imaging that is not present on early imaging would indicate neoplasm.

When the metastatic process is diffuse, virtually all of the radionuclide is concentrated in the skeleton, with little or no activity in the soft tissues or urinary tract. The resulting pattern, which is characterized by excellent bone detail, is frequently referred to as a *superscan*. A superscan may also be associated with metabolic bone diseases (e.g., Paget disease, renal osteodystrophy, or osteomalacia).

### Contraindications

- Patients who are pregnant, unless the benefits outweigh the risk of fetal damage
- Patients who are lactating

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for discussion of radiation exposure and risks.
- PT** Assure patients that they will not be exposed to large amounts of radioactivity because only tracer doses of the isotope are used.
- PT** Tell the patient that no fasting or sedation is required.
- PT** Inform the patient that the injection of the radionuclide may cause slight discomfort, nausea, or vomiting.

#### During

- Note the following procedural steps:
  1. The patient receives an IV injection of radionuclide, usually  $^{99m}\text{Tc}$ -MDP (technetium-99m-methylene diphosphonate) or  $^{99m}\text{Tc}$ -HDP 99m (hydroxymethane diphosphonate), into a peripheral vein in the arm.
  2. The patient is encouraged to drink several glasses of water between the time of radionuclide injection and the scanning. This facilitates renal clearance of the circulating tracer not picked up by the bone. The waiting period before scanning is approximately 2 to 3 hours.
  3. The patient is instructed to urinate.
  4. The patient is situated in the supine position on the scanning table in nuclear medicine.
  5. A radionuclide detector is placed over the patient's body and records the radiation emitted by the skeleton.
  6. This information is translated into a two- or three-dimensional view of the skeleton.

7. The patient may be repositioned in the prone and lateral positions during the test.
- Note that this scan is performed by a nuclear medicine technician in 30 to 60 minutes. It is interpreted by a physician trained in nuclear medicine imaging.
- PT** Inform patients in significant pain that lying on the hard scanning table can be uncomfortable.

### After

- Because only tracer doses of radionuclide are used, no radiation precautions need to be taken.
- PT** Assure the patient that the radioactive substance is usually excreted from the body within 6 to 24 hours.
- PT** Encourage the patient to drink fluids to aid in the excretion of the radioactive substance.

### Abnormal findings

Primary or metastatic tumor of the bone

Fracture

Degenerative arthritis

Rheumatoid arthritis

Osteomyelitis

Bone necrosis

Renal osteodystrophy

Paget disease of bone

### notes

**bone turnover markers** (BTMs, N-telopeptide [NTx], Bone collagen equivalents [BCEs], Osteocalcin [bone G1 a protein, BGP, osteocalc], Pyridinium [PYD] crosslinks, Bone-specific alkaline phosphatase [BSAP], Amino-terminal propeptide of type I procollagen [P1NP], C-telopeptide [CTx])

**Type of test** Blood; urine

### Normal findings

Results vary greatly with age.

#### N-telopeptide

Urine (nm BCE\*/mm creatinine)

Male: 21-83

Female, premenopausal: 17-94

Female, postmenopausal: 26-124

Serum (nm BCE\*)

Male: 5.4-24.2

Female: 6.2-19

(\*BCE = bone collagen equivalents)

#### C-telopeptide

Urine (ng/mL)

Adults:  $1.03 \pm 0.41$

Children:  $8.00 \pm 3.37$

Serum (pg/mL)

Female, premenopausal: 40-465

Female, postmenopausal: 104-1008

Male: 60-700

**Amino-terminal propeptide of type I procollagen, serum** (mcg/L)

Male: 22-105

Female, premenopausal: 19-101

Female, postmenopausal: 16-96

**Osteocalcin, serum** (ng/mL)

Adult (> 22 years)

Male: 5.8-14

Female: 3.1-14.4

**Pyridium, urine** (nm/mm)

Male: 10.3-33.6

Female: 15.3-33.6

**Bone-specific alkaline phosphatase, serum (mcg/L)**

Male: 6.5-20.1

Female, premenopausal: 4.5-16.9

Female, postmenopausal: 7-22.4

**Test explanation and related physiology**

With the increased use of bone density scans (see p. 157), osteoporosis can now be diagnosed and treated more easily. This has prompted an interest in biochemical markers of bone metabolism. Bone is continuously turned over: bone resorption by osteoclasts and bone formation by osteoblasts.

Bone mineral density studies (p. 157) are valuable tools in the identification of osteoporosis; however, they cannot recognize small changes in bone metabolism. Although bone density studies can be used to monitor the effectiveness of therapy, it takes years to detect measurable changes in bone density. Bone turnover markers (BTMs), however, can identify significant improvement in a few months after instituting successful therapy. Furthermore, the cost of bone density studies limits the feasibility of performing this test as frequently as may be required to monitor treatment.

Because the levels of BTMs vary according to the time of day and bone volume, these studies are not widely used or helpful in screening for detection of osteoporosis. Their use is in determining the effect of treatment by comparing BTMs with pretreatment levels. Levels will decline with the use of antiresorption drugs (e.g., estrogen, bisphosphonates, calcitonin, raloxifene). BTMs have been shown to be accurately predictive of early improvement in bone mineral density and antifracture treatment efficacy. BTMs are also useful in documenting treatment compliance.

*N*- and *C*-*telopeptides* (*NTx* and *CTx*) are protein fragments used in type I collagen that make up nearly 90% of the bone matrix. The *C* and *N* terminals of these proteins are crosslinked to provide tensile strength to the bone. When bone breaks down, *CTx* and *NTx* are released into the bloodstream and excreted in the urine. Serum levels of these fragments have been shown to correlate well with urine measurements normalized to creatinine. Measurements of these fragments show early response to anti-resorptive therapy (within 3-6 months) and are good indicators of bone resorption. Normal levels may vary with the method of testing.

## 170 bone turnover markers

*Amino-terminal propeptide of type I procollagen (PINP)*, like NTx, is directly proportional to the amount of new collagen produced by osteoblasts. Concentrations are increased in patients with various bone diseases that are characterized by increased osteoblastic activity. PINP is the most effective marker of bone formation and is particularly useful for monitoring bone formation therapies and antiresorptive therapies.

*Osteocalcin, or bone Gla protein (BGP)*, is a noncollagenous protein in the bone and is made by osteoblasts. It enters the circulation during bone resorption and bone formation; it is a good indicator of bone metabolism. Serum levels of BGP correlate with bone formation and destruction (turnover). Increased levels are associated with increased bone mineral density loss. BGP is a vitamin K-dependent protein. A reduced vitamin K intake is associated with reduced BGP levels. This probably explains the pathophysiology of vitamin K-dependent deficiency osteoporosis.

*Pyridinium (PYD) crosslinks* are formed during maturation of the type I collagen during bone formation. During bone resorption, these PYD crosslinks are released into the circulation.

*Bone-specific alkaline phosphatase (BSAP)* is an isoenzyme of alkaline phosphatase (see p. 29) and is found in the cell membrane of the osteoblast. It is, therefore, an indicator of the metabolic status of osteoblasts and bone formation.

These BTMs cannot indicate the risk of bone fracture nearly as well as a bone density measurement scan. These markers can be used to monitor the activity and treatment of Paget disease, hyperparathyroidism, and bone metastasis.

BTMs are normally high in children because of increased bone resorption associated with growth and remodeling of the ends of the long bones. The levels reach a peak at about age 14 years and then gradually decline to adult values. Because estrogen is a strong inhibitor of osteoclastic (bone resorption) activity, loss of bone density begins soon after menopause begins. Marker levels therefore rise after menopause. Most urinary assays are correlated with creatinine excretion for normalization.

### Interfering factors

- Measurements of these urinary markers can differ by as much as 30% in one person, even on the same day. Collecting double-voided specimens in the morning can minimize variability.
- Bodybuilding treatments, such as testosterone, can cause reduced levels of NTx.

## Procedure and patient care

### Blood

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: verify with laboratory
- Collect a venous blood sample in a red-top tube for NTx and/or a lavender- or green-top tube for osteocalcin. Check with the laboratory for guidelines with other markers.
- It is important to obtain baseline levels before instituting therapy.

### Urine

- See inside front cover for Routine Urine Testing.
- Preferably, obtain a double-voided specimen.
- Collect the first urine specimen 30 to 40 minutes before the time the specimen is needed.
- Discard this first specimen.
- Give the patient a glass of water to drink.
- At the requested time, obtain a second specimen.

## Abnormal findings

### ▲ Increased levels

Osteoporosis  
 Paget disease of bone  
 Advanced bone tumors  
 (primary or metastatic)  
 Acromegaly  
 Hyperparathyroidism  
 Hyperthyroidism

### ▼ Decreased levels

Hypoparathyroidism  
 Hypothyroidism  
 Cortisol therapy  
 Effective antiresorptive  
 therapy

## notes

**bone x-ray****Type of test** X-ray**Normal findings**

No evidence of fracture, tumor, infection, or congenital abnormalities

**Test explanation and related physiology**

X-ray images of the long bones are usually taken when the patient has complaints about a particular body area. Fractures or tumors are readily detectable by x-ray studies. In patients who have a severe or chronic infection overlying a bone (osteomyelitis), an x-ray image may detect the infection involving that bone. X-ray studies of the long bones also can detect joint destruction and bone spurring as a result of persistent arthritis. Growth patterns can be followed by serial x-ray studies of a long bone, usually the wrists and hands. Healing of a fracture can be documented and followed. X-ray images of the joints reveal the presence of joint effusions and soft tissue swelling as well. Calcifications in the soft tissue indicate chronic inflammatory changes of the nearby bursa or tendons. Soft tissue swelling also can be seen on these bone x-rays. Because the cartilage and tendons are not directly visualized, cartilage fractures, sprains, or ligamentous injuries cannot be seen.

At least two x-rays at 90-degree angles are required so that the bone region being studied can be visualized from two different angles (usually anterior to posterior and lateral). Some bone studies (e.g., skull, spine, hip) require oblique views to visualize all the parts that need to be seen.

**Interfering factors**

- Jewelry or clothing can obstruct radiographic visualization of part of the bone to be evaluated.
- Prior barium studies can diminish the full radiographic visualization of some of the bones surrounding the abdomen (e.g., spine and pelvis).

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- Handle carefully any injured parts of the patient's body.



- PT** Instruct the patient that he or she will need to keep the extremity still while the x-ray image is being taken. This can sometimes be difficult, especially when the patient has severe pain associated with a recent injury.
- Shield the patient's testes, ovaries, or pregnant abdomen to avoid exposure from scattered radiation.
- PT** Tell the patient that no fasting or sedation is required.

### During

- Note that, in the x-ray department, the patient is asked to place the involved extremity in several positions. An x-ray image is taken of each position.
  - Note that this test is routinely performed by a radiologic technologist within several minutes.
- PT** Tell the patient that no discomfort is associated with this test, except possibly from moving an injured extremity.

### After

- Administer an analgesic for relief of pain if indicated.

### Abnormal findings

Fractures

Congenital bone disorders (e.g., achondroplasia, dysplasia, dysostosis)

Tumors (osteogenic sarcoma, Paget disease, myeloma, or metastatic)

Infection or osteomyelitis

Osteoporosis or osteopenia

Joint destruction (arthritis)

Bone spurring

Abnormal growth pattern

Joint effusion

Foreign bodies

### notes

**brain scan** (Cisternogram, Cerebral blood flow, DaT scan)**Type of test** Nuclear scan**Normal findings**

No areas of altered radionuclide uptake within the brain

**Test explanation and related physiology**

The usefulness of a nuclear brain scan is narrow compared with computed tomography (CT), magnetic resonance imaging (MRI), and positron emission tomography (PET) scans of the brain. Primarily, a nuclear brain scan is used to indicate complete and irreversible cessation of brain function (brain death). This determination, when combined with appropriate clinical data, allows for cessation of medical therapy and opportunity for the harvest of potential donor organs. With brain death, there is complete absence of blood perfusion to the brain.

The brain scan can also be used to indicate cerebral vascular occlusion or stenosis. With the use of Diamox (acetazolamide), an accurate assessment of local cerebral blood flow can be determined. Diamox is a carbonic anhydrase inhibitor that results in the elevation of  $PCO_2$  in the bloodstream. Normally, this causes dilatation of the cerebral blood vessels. If asymmetric blood flow is noted after Diamox injection, cerebral vascular occlusion or stenosis can be suspected.

Brain scans are also used to investigate the ventricular system (cisternogram) of the central nervous system. Normal pressure hydrocephalus and ventricular shunt dysfunction can be identified and located.

**Contraindications**

- Patients who are pregnant unless the benefits outweigh the risk of fetal damage
- Patients who cannot cooperate during the testing

**Interfering factors**

- ✱ Many sedative drugs can affect brain nuclear imaging.
- ✱ Angiotensin-converting enzyme (ACE) inhibitors, vasoconstrictors, and vasodilators can alter blood flow distribution in nuclear brain imaging.

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.

- Administer blocking agents as ordered before scanning.
- Consider having a sedative ordered for agitated patients.

### During

- Note the following procedural steps:
  1. After administration of the radioisotope, the patient is placed in the supine position while planar and single-photon emission computed tomography (SPECT) images are obtained.
  2. When cerebral flow studies are performed, the counter is immediately placed over the head.
  3. The counts are anatomically recorded in timed sequence to follow the isotope during its first flow through the brain.
- Note that this study is performed by a technologist in the nuclear medicine department in approximately 35 to 45 minutes.

### After

- PT** Because only tracer doses of radioisotopes are used, inform the patient that no precautions need to be taken to prevent radioactive exposure to other personnel or family present.
- PT** Encourage the patient to drink fluids to aid in the excretion of the isotope from the body.

### Abnormal findings

Cerebral death

Cerebral vascular stenosis/occlusion

Cerebral neoplasm

CSF leakage

Hydrocephalus

### notes

**breast cancer tumor analysis** (DNA ploidy status, S-phase fraction, Cathepsin D, HER-2 [c erbB2, neu] protein, p53 protein, Ki67 protein)

**Type of test** Microscopic examination

### Normal findings

#### DNA ploidy

Aneuploid is unfavorable.

Diploid is favorable.

#### S-phase fraction

> 5.5% is unfavorable.

< 5.5% is favorable.

#### HER-2 protein

IHC method: 0 to 1 +

FISH method: < 2 copies/cell

Oncotype DX method: < 10.7 units

#### Cathepsin D

> 10% is unfavorable.

< 10% is favorable.

#### p53 protein

> 10% is unfavorable.

< 10% is favorable.

#### Ki67 protein

> 20% is unfavorable.

10% to 20% is borderline.

< 20% is favorable.

### Test explanation and related physiology

The most important predictor of recurrent breast cancer is stage of disease, including lymph node status. Patients with positive lymph node metastasis are more likely to develop recurrence. However, nearly 30% of the patients whose tumors have been completely removed and who have no evidence of lymph node metastasis will also develop recurrence. Conventional predictors such as tumor size, grade, histologic type, and hormone receptors (see p. 389) are able to identify some of the patients who are at increased risk for recurrence. However, it is important to accurately predict the patients who are destined for recurrence so that they can be selected for systemic therapy; patients who will not have a recurrence can be spared the morbidity of a treatment that is not needed.

## DNA ploidy status and S-phase fraction

These are measurements of the rapidity with which the cells in a breast cancer grow.

## Cathepsin D

This protein catabolic enzyme was found to be absent in resting breast tissue but significantly elevated in malignant tissue. This protein exists on tumor cell membrane and is correlated with worse clinical outcomes. The exact cutoff point between a favorable prognosis and unfavorable prognosis has yet to be standardized.

## HER-2 (c erbB2, neu) protein

*HER-2/neu*, which stands for *human epidermal growth factor receptor 2*, is a protein associated with worse clinical outcomes. The *HER-2/neu* oncogene encodes a transmembrane tyrosine kinase receptor with extensive similarity to other epidermal growth factor receptors. It is normally involved in the pathways leading to cell growth and survival. Approximately 15% to 20% of breast cancers have an amplification of the *HER-2/neu* gene or overexpression of its protein product.

The *HER-2* gene can act as a target for an antineoplastic monoclonal antibody drug such as trastuzumab (Herceptin).

## p53 protein

The *p53* gene is a tumor suppressor gene that is overexpressed in more aggressive breast cancer cells. Mutation of the gene causes overexpression and a buildup of mutant proteins on the surface of the cancer cells.

## Ki67 protein

The *Ki67* gene encodes the synthesis for the Ki67 protein that is associated with worse clinical outcomes.

## Interfering factors

- Delay in tissue fixation may cause deterioration of marker proteins and may produce lower values.
- ✱ Preoperative use of some chemotherapy agents may cause *decreased* levels of some marker proteins.

## Procedure and patient care

### Before

- PT** Inform the patient that an examination for these tumor predictor markers may be performed on their breast cancer tissue.
- PT** Provide psychological and emotional support to the breast cancer patient.

## 178 breast cancer tumor analysis

### During

- The surgeon obtains tumor tissue.
- This tissue should be placed on ice or in formalin.
- Part of the tissue is used for routine histology. A portion of the paraffin block is sent to a reference laboratory.

### After

- ▶ Explain to the patient that results are usually available in 1 week.

### Abnormal findings

Unfavorable test results indicating a risk of cancer reoccurrence

notes

## breast ductal lavage

**Type of test** Fluid analysis

### Normal findings

No atypical cells in the effluent

### Possible critical values

Cancer cells in the effluent

### Test explanation and related physiology

The theory behind ductal lavage is that by washing out exfoliated cells from a few breast ducts, the risk of breast cancer developing in the near future can be assessed. If atypical cells are obtained, the risk of breast cancer developing in the next decade may be as high as 4 to 10 times normal. Once that risk is identified, the patient may choose to attempt to alter that risk by using chemopreventative medications (e.g., selective estrogen receptor modulators) or surgery.

Its use is limited to women who have been found to be at a statistically higher risk for breast cancer by *Gail or Claus breast cancer risk models*. These statistical models are based on age at menarche, age at first pregnancy, prior breast surgery, family history, and history of atypical changes in previous breast biopsies.

### Contraindications

- Patients with prior breast cancer surgery because their risks are already known to be high

### Potential complications

- Infection

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. Often these women have already received extensive counseling regarding their risk of breast cancer.
- Be sure that the breast examination and mammogram are normal.
  - Apply a topical anesthetic to the nipple area about a half hour before the test.

## 180 breast ductal lavage

### During

- Note the following procedural steps:
  1. A suction apparatus is applied to the nipple area. Ducts that reveal fluid with the suction are then chosen for cannulation.
  2. A tiny catheter is gently placed into the nipple duct, and the duct is lavaged with 1 to 3 mL of saline.
  3. The effluent is collected in a small tube and sent for cytology.
  4. The procedure is then repeated for the other ducts that produced fluid with nipple suction.
- This procedure is performed by a surgeon in the office. There is minimal to moderate discomfort associated with the nipple suction, duct cannulation, and lavage.

### After

- PT** Inform the patient of the possibility of mild breast discomfort.
- Arrange for follow-up to discuss test results.
- PT** Provide counseling if results indicate atypical or malignant cells.

### Abnormal findings

Atypical cells

Ductal cancer cells

### notes



## breast ultrasonography (Ultrasound mammography, Breast sonogram)

B

**Type of test** Ultrasound

### Normal findings

No evidence of cyst or tumor

### Test explanation and related physiology

Ultrasound examination of the breast is diagnostically performed to determine whether a mammographic abnormality or a palpable lump is a cyst (fluid-filled) or a solid tumor (benign or malignant). It is also used in screening for breast cancer in women whose breasts are dense on mammography.

Ultrasonography of the breast is used to:

- Differentiate cystic from solid breast lesions.
- Identify masses in women with breast tissue too dense for accurate mammography.
- Monitor a cyst to determine whether it enlarges or disappears.
- Measure the size of a tumor.
- Evaluate the axilla in women who are newly diagnosed with breast cancer.

Ultrasonography is also useful for examination of symptomatic breasts in women in whom the radiation of mammography is potentially harmful. These include:

- Pregnant women. Radiation may be harmful to the fetus.
- Women younger than age 25 years, who may be at greater oncologic risk from the radiation of mammography.
- Women who refuse mammography because of unreasonable fear of diagnostic radiation.

With high-quality diagnostic ultrasonography, the characteristics of an abnormality can be evaluated and a reasonable prediction can be made whether it is malignant. Diagnostic accuracy is improved when breast ultrasonography is combined with mammography (see p. 606).

Ultrasound can be used to locate and accurately direct percutaneous biopsy probes to a nonpalpable breast abnormality for biopsy or aspiration. Ultrasound is painless, harmless, and without any radiation effects on the breast tissue.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

## 182 breast ultrasonography

- PT** Assure the patient that no discomfort is associated with this study.
- PT** Inform the patient that no fasting or sedation is required before the tests. Instruct the patient not to apply any lotions or powders to the breasts on the examination day.

### During

- The patient is placed in the supine position, and the transducer is directly applied to the breast using contact gel to improve sound transmission.
- Note that this test is performed by an ultrasound technician in approximately 15 minutes.
- There is no discomfort associated with this procedure.

### After

- After the test is completed, the breasts are dried, and the conductive paste is removed.

### Abnormal findings

Cyst

Hematoma

Cancer

Fibroadenoma

Fibrocystic disease

Abscess

### notes

## bronchoscopy

**Type of test** Endoscopy

### Normal findings

Normal larynx, trachea, bronchi, and alveoli

### Test explanation and related physiology

Bronchoscopy permits endoscopic visualization of the larynx, trachea, and bronchi by either a flexible fiberoptic bronchoscope or a rigid bronchoscope. There are many diagnostic and therapeutic uses for bronchoscopy.

*Diagnostic* uses of bronchoscopy include:

- Direct visualization of the tracheobronchial tree for abnormalities (e.g., tumors, inflammation, strictures)
- Biopsy of tissue from observed lesions
- Aspiration of deep sputum for culture, sensitivity, and cytology determinations
- Direct visualization of the larynx for identification of vocal cord paralysis if present

*Therapeutic* uses of bronchoscopy include:

- Aspiration of retained secretions in patients with airway obstruction or postoperative atelectasis
- Control of bleeding within the bronchus
- Removal of foreign bodies that have been aspirated
- Brachytherapy, which is endobronchial radiation therapy using an iridium wire placed via the bronchoscope
- Palliative laser obliteration of bronchial neoplastic obstruction
- Access for ultrasound

The *flexible fiberoptic bronchoscope* has accessory lumens through which cable-activated instruments can be used for removing biopsy specimens of pathologic lesions. Also, the collection of bronchial washings (obtained by flushing the airways with saline solution), pulmonary toilet, and the instillation of anesthetic agents can be carried out through these extra lumens. Double-sheathed, plugged-protected brushes also can be passed through this accessory lumen. Specimens for cytology and bacteriology can be obtained with these brushes.

Needles can be placed through the scope to obtain biopsies from tissue immediately adjacent to the bronchi. Transbronchial needle aspiration can be directed by the use of *endobronchial ultrasound*. This technique is particularly helpful in staging lung cancers and identifying sarcoidosis, lymphomas, and infections.

## 184 bronchoscopy

Laser therapy to burn out endotracheal lesions can now be performed through the bronchoscope.

*Laryngoscopy* is often performed through a short bronchoscope to allow inspection of the larynx and paralaryngeal structures. This is most commonly performed by an ENT (ear, nose, and throat) surgeon. Cancers, polyps, inflammation, and infections of these structures can be identified. The vocal cord motion can be evaluated also. Anesthesiologists use laryngoscopy to visualize the vocal cord structures on patients who are difficult to intubate for general anesthesia.

### Contraindications

- Patients with hypercapnia and severe shortness of breath who cannot tolerate interruption of high-flow oxygen
- Severe tracheal stenosis may make it difficult to pass the scope.

### Potential complications

- Fever
- Hypoxemia
- Laryngospasm
- Bronchospasm
- Pneumothorax
- Aspiration
- Hemorrhage (after biopsy)


### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. Allay any fears and allow the patient to verbalize any concerns.
- Obtain informed consent for this procedure.
- Keep the patient NPO for 4 to 8 hours before the test to reduce the risk of aspiration.
- PT** Instruct the patient to perform good mouth care to minimize the risk of introducing bacteria into the lungs.
- Remove and safely store the patient's dentures, glasses, or contacts before administering the preprocedure medications.
- Administer the preprocedure medications as ordered.
- PT** Reassure the patient that he or she will be able to breathe during this procedure.
- PT** Instruct the patient not to swallow the local anesthetic sprayed into the throat. Provide a basin for expectoration of the lidocaine.
- Have emergency resuscitation equipment available.

## During

- Note the following procedural steps for *fiberoptic bronchoscopy*:
  1. This test is performed by a pulmonary specialist or a surgeon at the bedside or in an appropriately equipped room.
  2. The patient's nasopharynx and oropharynx are anesthetized topically with lidocaine spray before insertion of the bronchoscope.
  3. The patient is placed in a sitting or supine position, and the tube is inserted through the nose or mouth and into the pharynx (Figure 7).
  4. After the tube is passed into the larynx and through the glottis, more lidocaine is sprayed into the trachea to prevent the cough reflex.

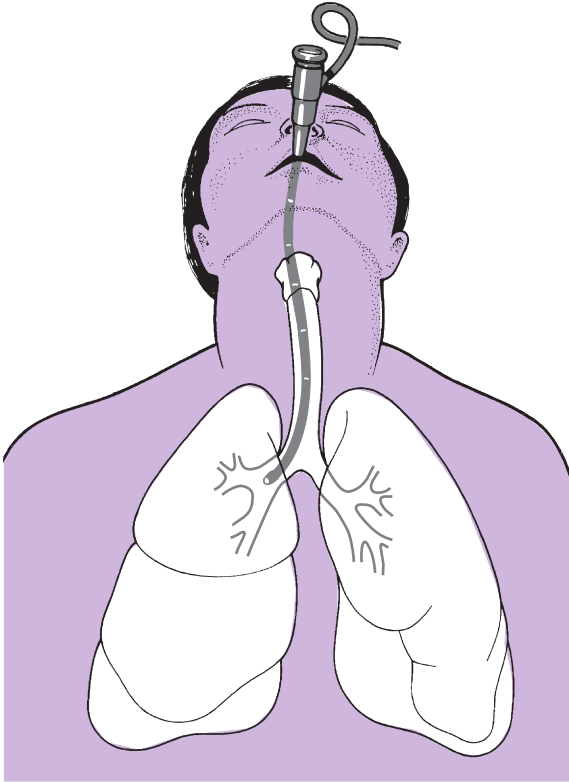

**FIG. 7** Bronchoscopy. A bronchoscope is inserted through the trachea and into the bronchus.

## 186 bronchoscopy

5. The tube is passed farther, well into the trachea, bronchi, and first- and second-generation bronchioles, for systematic examination of the bronchial tree.
  6. Biopsy specimens and washings are taken if pathology is suspected.
- Note that this procedure is performed by a physician in approximately 30 to 45 minutes.
- PT** Tell the patient that because of sedation, no discomfort is usually felt.

### After

- PT** Instruct the patient not to eat or drink anything until the tracheobronchial anesthesia has worn off and the gag reflex has returned, usually in approximately 2 hours.
- Observe the patient's sputum for hemorrhage if biopsy specimens were removed. A small amount of blood streaking may be expected and is normal for several hours after the procedure. Large amounts of bleeding can cause a chemical pneumonitis.
  - Observe the patient closely for evidence of impaired respiration or laryngospasm. The vocal cords may go into spasms after intubation.
- PT** Inform the patient that postbronchoscopy fever often develops within the first 24 hours. High, persistent fever should be reported immediately.
- If a tumor is suspected, collect a postbronchoscopy sputum sample for a cytology determination.
- PT** Inform the patient that warm saline gargles and lozenges may be helpful if a sore throat develops.
- PT** Inform the patient that biopsy or culture reports will be available in 2 to 7 days.

### Abnormal findings

Inflammation  
Strictures  
Tuberculosis  
Cancer  
Hemorrhage  
Foreign body  
Abscess  
Infection

### notes

**calcitonin** (Human calcitonin [HCT], Thyrocalcitonin)**Type of test** Blood**Normal findings**

Basal (plasma)

Males:  $\leq 19$  pg/mL or  $\leq 19$  ng/L (SI units)Females:  $\leq 14$  pg/mL or  $\leq 14$  ng/L (SI units)

Calcium infusion (2.4 mg/kg)

Males:  $\leq 190$  pg/mL or  $\leq 190$  ng/LFemales:  $\leq 130$  pg/mL or  $\leq 130$  ng/L

Pentagastrin injection (0.5 mcg/kg)

Males:  $\leq 110$  pg/mL or  $\leq 110$  ng/LFemales:  $\leq 30$  pg/mL or  $\leq 30$  ng/L**Test explanation and related physiology**

Calcitonin is a hormone secreted by the parafollicular or C-cells of the thyroid gland. Its secretion is stimulated by elevated serum calcium levels. The purpose of calcitonin is to contribute to calcium homeostasis.

This test is usually used in the evaluation of patients with or suspected to have medullary carcinoma of the thyroid (see p. 195). This is a cancer of the thyroid with a familial tendency; if it is found late, it has a poor prognosis. Routine screening for elevated calcitonin levels can detect medullary cancer early and can improve chances for cure. C-cell hyperplasia, a benign calcitonin-producing disease that also has a familial tendency, is also associated with elevated calcitonin levels.

Equivocal elevations in calcitonin levels should be followed with further provocative testing using pentagastrin or calcium to stimulate calcitonin secretion. *Pentagastrin stimulation* involves an intravenous (IV) infusion with blood samples drawn before the injection and at 90 seconds, 2 minutes, and 5 minutes after the infusion. The *calcium infusion test* can be performed in a variety of ways but is most commonly administered with baseline and 5- and 10-minute postinfusion blood levels.

Elevated levels of calcitonin also may be seen in people with cancer of the lung, breast, or pancreas. This is probably a form of paraneoplastic syndrome in which there is an ectopic production of calcitonin by the nonthyroid cancer cells.

**Interfering factors**

- Levels are often elevated in pregnancy and in newborns.

## 188 calcitonin

- ❖ Drugs that may cause *increased* levels include calcium, cholecystokinin, epinephrine, glucagon, pentagastrin, and oral contraceptives.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: yes
  - Blood tube commonly used: green or red
  - Collect a venous sample of blood in a heparinized green-top tube or a chilled red-top tube according to the laboratory's protocol.
  - The specimen should be placed on ice immediately. The blood may be frozen and sent to a reference laboratory.
- PT** Tell the patient that results may not be available for several days.

### Abnormal findings

#### ▲ Increased levels

Medullary carcinoma of the thyroid  
C-cell hyperplasia  
Oat cell carcinoma of the lung  
Breast carcinoma  
Pancreatic cancer  
Primary hyperparathyroidism  
Secondary hyperparathyroidism because of chronic renal failure  
Pernicious anemia  
Zollinger–Ellison syndrome  
Alcoholic cirrhosis  
Thyroiditis

### notes



**calcium** (Total/ionized calcium, Ca, Serum calcium)**Type of test** Blood; urine**Normal findings**

Age	mg/dL	mmol/L
<b>Total calcium</b>		
< 10 days	7.6-10.4	1.9-2.6
Umbilical	9-11.5	2.25-2.88
10 days-2 years	9-10.6	2.3-2.65
Child	8.8-10.8	2.2-2.7
Adult*	9-10.5	2.25-2.62
<b>Ionized calcium</b>		
Newborn	4.2-5.58	1.05-1.37
2 months-18 years	4.8-5.52	1.2-1.38
Adult	4.5-5.6	1.05-1.3

\* In elderly individuals, values tend to decrease.

**Possible critical values**

Total calcium: &lt; 6 or &gt; 13 mg/dL

Ionized calcium: &lt; 2.2 or &gt; 7 mg/dL

**Test explanation and related physiology**

The serum calcium test is used to evaluate parathyroid function and calcium metabolism by directly measuring the total amount of calcium in the blood. Determination of serum calcium is used to monitor patients with renal failure, renal transplantation, hyperparathyroidism, and various malignancies. It is also used to monitor calcium levels during and after large-volume blood transfusions.

About half the total calcium in the blood exists in its free (ionized) form, and about half exists in its protein-bound form (mostly with albumin). The serum calcium level is a measure of both. As a result, when the serum albumin level is low (as in malnourished patients), the serum calcium level will also be low and vice versa. As a rule of thumb, the total serum calcium level decreases by approximately 0.8 mg for every 1-g decrease in the serum albumin level. Serum albumin should be measured with serum calcium.

When the serum calcium level is elevated on at least three separate determinations, the patient is said to have *hypercalcemia*.

The most common cause of hypercalcemia is hyperparathyroidism. Parathyroid hormone (PTH; see p. 675) causes elevated calcium levels by increasing gastrointestinal (GI) absorption, decreasing urinary excretion, and increasing bone resorption. Malignancy, the second most common cause of hypercalcemia, can cause elevated calcium levels in two main ways. First, tumor metastasis (myeloma, lung, breast, renal cell) to the bone can destroy the bone, causing resorption and pushing calcium into the blood. Second, the cancer (lung, breast, renal cell) can produce a PTH-like substance that drives the serum calcium up (ectopic PTH). Excess vitamin D ingestion can increase serum calcium by increasing renal and GI absorption. Granulomatous infections, such as sarcoidosis and tuberculosis, are associated with hypercalcemia.

Hypocalcemia occurs in patients with hypoalbuminemia. The most common causes of hypoalbuminemia are malnutrition (especially in alcoholics) and large-volume IV infusions. Large blood transfusions are associated with low serum calcium levels because the citrate additives used in banked blood for anticoagulation bind the free calcium in the recipient's bloodstream. Intestinal malabsorption, renal failure, rhabdomyolysis, alkalosis, and acute pancreatitis (caused by saponification of fat) are also known to be associated with low serum calcium levels.

*Urinary calcium* can also be measured. Excretion of calcium in the urine is increased in all patients with hypercalcemia. Urinary calcium levels are decreased in patients with hypocalcemia. The test is helpful in determining the cause of recurrent nephrolithiasis.

### Interfering factors

- Vitamin D intoxication may cause increased calcium levels.
  - Excessive ingestion of milk may cause increased levels.
  - Serum pH can affect calcium values. A decrease in pH causes increased calcium levels.
  - Prolonged tourniquet time will lower pH and falsely increase calcium levels.
  - There is normally a small diurnal variation in calcium, with peak levels occurring around 9 PM.
  - Hypoalbuminemia is artifactually associated with decreased levels of total calcium.
- ▼ Drugs that may cause *increased* serum levels include alkaline antacids, androgens, calcium salts, ergocalciferol, hydralazine, lithium, PTH, thiazide diuretics, thyroid hormone, and vitamin D.

- Drugs that may cause *decreased* serum levels include acetazolamide, albuterol, anticonvulsants, asparaginase, aspirin, calcitonin, cisplatin, corticosteroids, diuretics, estrogens, heparin, laxatives, loop diuretics, magnesium salts, and oral contraceptives.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: verify with laboratory
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels (hypercalcemia)

Hyperparathyroidism  
 Nonparathyroid PTH-producing tumor (e.g., lung or renal carcinoma)  
 Metastatic tumor to the bone  
 Paget disease of bone  
 Prolonged immobilization  
 Milk-alkali syndrome  
 Vitamin D intoxication  
 Lymphoma  
 Granulomatous infections (e.g., sarcoidosis and tuberculosis)  
 Addison disease  
 Acromegaly  
 Hyperthyroidism

#### ▼ Decreased levels (hypocalcemia)

Hypoparathyroidism  
 Renal failure  
 Hyperphosphatemia secondary to renal failure  
 Rickets  
 Vitamin D deficiency  
 Osteomalacia  
 Malabsorption  
 Pancreatitis  
 Fat embolism  
 Alkalosis  
 Hypoalbuminemia

### notes

**caloric study** (Oculovestibular reflex study)**Type of test** Electrodiagnostic**Normal findings**

Nystagmus with irrigation

**Test explanation and related physiology**

Caloric studies are used to evaluate the vestibular portion of the eighth cranial nerve (CN VIII) by irrigating the external auditory canal with hot or cold water. This is considered a part of a complete neurologic examination. Normally stimulation with cold water causes rotary nystagmus (involuntary rapid eye movement) away from the ear being irrigated; hot water induces nystagmus toward the side of the ear being irrigated. If the labyrinth is diseased or CN VIII is not functioning (e.g., from tumor compression), no nystagmus is induced. This study aids in the differential diagnosis of abnormalities that may occur in the vestibular system, brainstem, or cerebellum.

**Contraindications**

- Patients with a perforated eardrum
- Cold air may be substituted for the fluid, although this method is much less reliable.
- Patients with an acute disease of the labyrinth (e.g., Ménière syndrome)
- The test can be performed when the acute attack subsides.

**Interfering factors**

- ✚ Drugs such as sedatives and antivertigo agents can alter test results.

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to avoid solid foods before the test to reduce the incidence of vomiting.

**During**

- Although the exact procedures for caloric studies vary, note the following steps in a typical test:
  1. Before the test, the patient is examined for the presence of nystagmus, postural deviation (Romberg sign), and past-pointing. This examination provides the baseline values for comparison during the test.

2. The ear canal should be examined and cleaned before testing to ensure that the water will freely flow to the middle ear area.
  3. The ear on the suspected side is irrigated first because the patient's response may be minimal.
  4. After an emesis basin is placed under the ear, the irrigation solution is directed into the external auditory canal until the patient complains of nausea and dizziness or nystagmus is seen. Usually this occurs in 20 to 30 seconds.
  5. If after 3 minutes no symptoms occur, the irrigation is stopped.
  6. The patient is tested again for nystagmus, past-pointing, and Romberg sign.
  7. After approximately 5 minutes, the procedure is repeated on the other side.
- Note that this procedure is usually performed by a physician or technician in approximately 15 minutes.
- PT** Tell the patient that he or she will probably experience nausea and dizziness during the test.

### After

- Usually place the patient on bed rest for approximately 30 to 60 minutes until nausea or vomiting subsides.
- Ensure patient safety related to dizziness.

### Abnormal findings

Brainstem inflammation, infarction, or tumor

Cerebellar inflammation, infarction, or tumor

Vestibular or cochlear inflammation or tumor

Acoustic neuroma

CN VIII neuritis or neuropathy

### notes

**cancer tumor markers** (Tumor markers [TMs], Tumor-associated markers)

**Type of test** Blood

**Normal findings**

Normal values vary per laboratory

**Test explanation and related physiology**

Tumor markers are produced by cancer cells or by other cells of the body in response to cancer. Most TMs are made by normal cells as well as by cancer cells; however, they are overproduced by cancers. These substances can be found in the blood, urine, stool, tumor tissue, or other tissues or bodily fluids. Most TMs are proteins. However, changes in tumor cell RNA and DNA are also used as TMs. Many different TMs have been characterized. Some are associated with only one type of cancer, but others are associated with two or more cancer types. Sometimes noncancerous conditions can cause elevated TMs, but they are usually less elevated than in cancer.

Tumor markers are generally not used for cancer screening because they are not sensitive or specific enough. Neither are they used to diagnose cancers. Prostate-specific antigen (PSA), CA-125, and a few other TMs are being used very carefully in screening, but, alone, the cost-effectiveness is questionable.

Tumor marker levels may be measured before treatment to help doctors plan the appropriate therapy. In some types of cancer, the level of a tumor marker reflects the stage of the disease and/or the patient's prognosis. Pretreatment TM scan be used to establish a baseline level against which posttreatment TMs can be compared in order to determine the effectiveness of therapy. A decrease in the level of a TM or a return to the TM's normal level may indicate that the cancer is responding to treatment, but no change or an increase may indicate that the cancer is not responding.

Tumor markers may also be measured during cancer follow-up evaluations to check for recurrent disease. There are many TMs presently being used in cancer care. [Table 6](#) lists the most commonly used markers.

**Interfering factors**

- Other benign and malignant diseases associated with elevations

**TABLE 6** Tumor markers and their associated cancers

Tumor marker	Associated cancer
<i>ALK</i> gene	Non-small cell lung cancer Lymphoma
Alpha-fetoprotein (AFP)	Liver Germ cell tumors
BCR-ABL fusion gene (Philadelphia chromosome)	Chronic myeloid leukemia Acute lymphoblastic leukemia Acute myelogenous leukemia
Beta 2 microglobulin	Liver Germ cell tumors
BRAF V600 mutations	Cutaneous melanoma Colorectal
C-kit/CD117	Gastrointestinal stromal tumor Melanoma
CA15-3/CA27.29 CA19-9	Breast Pancreas Biliary Stomach
CA-125	Ovary
Calcitonin	Medullary thyroid carcinoma
Carcinoembryonic antigen (CEA)	Colon Other gastrointestinal tumors Breast
CD20	Non-Hodgkin lymphoma
Chromogranin A	Neuroendocrine
Chromosomes 3, 17, and 9p21	Bladder
Cytokeratin fragment 21-1	Lung
Des-gamma-carboxy prothrombin	Liver
Fibrin/fibrinogen	Bladder
HE4	Ovary
Human chorionic gonadotropin (HCG)	Choriocarcinoma Germ cell tumors
Immunoglobulins	Multiple myeloma Waldenström macroglobulinemia
Inhibin A	Germ cell tumors of ovary

*Continued*

**TABLE 6** Tumor markers and their associated cancers—cont'd

Tumor marker	Associated cancer
KRAS	Colorectal Non-small cell lung cancer
Lactate dehydrogenase (LDH)	Germ cell tumors Leukemia Melanoma Brain
Neuron-specific enolase (NSE)	Small cell lung cancer Neuroblastoma
Nuclear matrix protein 22	Bladder
Plasminogen activator inhibitor (PAI-1)	Breast
Programmed death ligand 1 (PD-L1)	Non-small cell lung cancer
Prostate-specific antigen (PSA)	Prostate
Squamous cell carcinoma (SCC) antigen	Squamous cell carcinoma of the cervix, oral cavity, esophagus, lung, anal canal, and skin
Thyroglobulin	Thyroid
Urokinase plasminogen activator (uPA)	Breast

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: varies by test and laboratory
- The blood sample may be sent to a central diagnostic laboratory. The results may not be available for 7 to 10 days.

### Abnormal findings

#### ▲ Increased levels

- Cancer
- Metastatic cancer
- Benign disease



**carbon dioxide content** (CO<sub>2</sub> content, CO<sub>2</sub> combining power)**Type of test** Blood**Normal findings**

Adult/elderly: 23-30 mEq/L or 23-30 mmol/L (SI units)

Child: 20-28 mEq/L

Infant: 20-28 mEq/L

Newborn: 13-22 mEq/L

**Possible critical values**

&lt; 6 mEq/L

**Test explanation and related physiology**

The CO<sub>2</sub> content is a measure of CO<sub>2</sub> in the blood. In the peripheral venous blood, this assists in evaluation of the pH status of the patient and in evaluation of electrolytes. The serum CO<sub>2</sub> test is usually included with other assessments of electrolytes. It is usually done with a multiphasic testing machine that also measures sodium, potassium, chloride, blood urea nitrogen (BUN), and creatinine.

It is important not to get this test confused with Pco<sub>2</sub>. This CO<sub>2</sub> content measures the H<sub>2</sub>CO<sub>3</sub>, the dissolved CO<sub>2</sub>, and the bicarbonate ion (HCO<sub>3</sub><sup>-</sup>) that exists in the serum. Because the amounts of H<sub>2</sub>CO<sub>3</sub> and dissolved CO<sub>2</sub> in the blood are so small, CO<sub>2</sub> content is an indirect measure of the HCO<sub>3</sub><sup>-</sup> anion. The HCO<sub>3</sub><sup>-</sup> anion is second in importance to the chloride ion in electrical neutrality (negative charge) of extracellular and intracellular fluid; its major role is in acid–base balance.

Levels of HCO<sub>3</sub><sup>-</sup> are regulated by the kidneys. Increases cause alkalosis, and decreases cause acidosis. See further discussion of this test as it is performed on arterial blood (see p. 106). When CO<sub>2</sub> content is measured in the laboratory with other serum electrolytes, air affects the specimen, and the CO<sub>2</sub> partial pressure can be altered. Therefore venous blood specimens are not very accurate for true CO<sub>2</sub> content or HCO<sub>3</sub><sup>-</sup> determination. This test is used mostly as a rough guide to the patient's acid–base balance.

**Interfering factors**

- Underfilling the tube with blood allows CO<sub>2</sub> to escape from the serum specimen and may significantly decrease HCO<sub>3</sub><sup>-</sup> values.
- ✚ Drugs that may cause *increased* serum CO<sub>2</sub> and HCO<sub>3</sub><sup>-</sup> levels include aldosterone, barbiturates, bicarbonates, ethacrynic acid, hydrocortisone, loop diuretics, mercurial diuretics, and steroids.

## 198 carbon dioxide content

- ▼ Drugs that may cause *decreased* levels include methicillin, nitrofurantoin (Furadantin), paraldehyde, phenformin hydrochloride, tetracycline, thiazide diuretics, and triamterene.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or green

### Abnormal findings

#### ▲ Increased levels

Severe diarrhea  
Starvation  
Severe vomiting  
Aldosteronism  
Emphysema  
Metabolic alkalosis  
Gastric suction

#### ▼ Decreased levels

Renal failure  
Salicylate toxicity  
Diabetic ketoacidosis  
Metabolic acidosis  
Shock  
Starvation

### notes

**carboxyhemoglobin (COHb, Carbon monoxide)****Type of test** Blood**Normal findings**

Saturation of hemoglobin

Nonsmoker: &lt; 3%

Smoker:  $\leq$  12%Newborn:  $\geq$  12%**Possible critical values**

&gt; 20%

**Test explanation and related physiology**

This test is used to detect carbon monoxide poisoning. It measures the amount of serum COHb, which is formed by the combination of carbon monoxide (CO) and hemoglobin (Hb). CO combines with Hb 200 times more readily than oxygen (O<sub>2</sub>) can combine with Hb; thus fewer Hb bonds are available to combine with O<sub>2</sub>. Furthermore, when CO occupies the O<sub>2</sub> binding sites, Hb is changed to bind the remaining O<sub>2</sub> more tightly. This greater affinity of CO for Hb and this change in O<sub>2</sub> binding strength do not allow O<sub>2</sub> to pass readily from RBCs to tissue. Less O<sub>2</sub> is therefore available for tissue cell respiration. This results in hypoxemia.

Carbon monoxide poisoning is documented by Hb analysis for COHb. A specimen should be drawn as soon as possible after exposure because CO is rapidly cleared from Hb by breathing normal air. O<sub>2</sub> saturation studies and oximetry are inaccurate in CO-exposed patients because they measure all forms of oxygen-saturated Hb, including COHb. In these circumstances, the patient's oximetry will be good, yet the patient will be hypoxicemic.

This test can also be used to evaluate patients with complaints of headache, irritability, nausea, vomiting, and vertigo, who unknowingly may have been exposed to CO. Its greatest use, however, is in patients exposed to smoke inhalation, exhaust fumes, and fires. Other sources of CO include tobacco smoke, petroleum and natural gas fuel fumes, automobile exhaust, unvented natural gas heaters, and defective gas stoves.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender or green
- Obtain the patient's history for possible sources of CO.

## 200 carboxyhemoglobin

- Assess the patient for signs and symptoms of mild CO toxicity (e.g., headache, weakness, dizziness, malaise, dyspnea) and moderate to severe CO toxicity (e.g., severe headache, bright red mucous membranes, cherry-red blood). Maintain patient safety precautions if confusion is present.
  - Treat the patient as indicated by the physician. Usually the patient receives high concentrations of O<sub>2</sub> to displace the COHb.
- PT** Encourage respirations to allow the patient to clear CO from the Hb.

### Abnormal findings

Carbon monoxide poisoning

notes

**cardiac catheterization** (Coronary angiography, Ventriculography)

**Type of test** X-ray with contrast dye

**Normal findings**

Normal heart-muscle motion, normal coronary arteries, normal great vessels, and normal intracardiac pressures and volumes

**Test explanation and related physiology**

Cardiac catheterization is used to visualize the heart chambers, arteries, and great vessels. It is used most often to evaluate patients with chest pain. Patients with positive stress test results are also studied to locate the region of coronary occlusion. This test is also used to determine the effects of valvular heart disease. Right heart catheterization is performed to calculate cardiac output and to measure right heart pressures. Right heart catheterization is also used to identify pulmonary emboli (see pulmonary angiography, p. 757).

For cardiac catheterization, a catheter is passed into the heart through a peripheral vein or artery, depending on whether catheterization of the right or left side of the heart is being performed. Pressures are recorded through the catheter, and radiographic dyes are injected. Cardiac output and other measures of cardiac functions can be determined. Cardiac catheterization is indicated for the following reasons:

- To identify, locate, and quantitate the severity of atherosclerotic, occlusive coronary artery disease
- To evaluate the severity of acquired and congenital cardiac valvular or septal defects
- To detect congenital cardiac abnormalities, such as transposition of great vessels, patent ductus arteriosus, and anomalous venous return to the heart
- To evaluate the success of previous cardiac surgery or balloon angioplasty
- To evaluate cardiac muscle function
- To identify and quantify ventricular aneurysms
- To detect disease of the great vessels, such as atherosclerotic occlusion or aneurysms within the aortic arch
- To evaluate and treat patients with acute myocardial infarction
- To insert a catheter to monitor right-sided heart pressures, such as pulmonary artery and pulmonary wedge pressures (Table 7 provides pressures and volumes used in cardiac monitoring.)

**TABLE 7** Pressures and volumes used in cardiac monitoring

Pressures/volumes	Description	Normal values
<b>Pressures</b>		
Routine blood pressure	Routine brachial artery pressure	90-140/60-90 mm Hg
Systolic left ventricular pressure	Peak pressure in the left ventricle during systole	90-140 mm Hg
End-diastolic left ventricular pressure	Pressure in the left ventricle at the end of diastole	4-12 mm Hg
Central venous pressure	Pressure in the superior vena cava	2-14 cm H <sub>2</sub> O
Pulmonary wedge pressure	Pressure in the pulmonary venules, an indirect measurement of left atrial pressure and left ventricular end-diastolic pressure	Left atrial: 6-15 mm Hg
Pulmonary artery pressure	Pressure in the pulmonary artery	15-28/5-16 mm Hg
Aortic artery pressure	Same as routine blood pressure	
<b>Volumes</b>		
End-diastolic volume (EDV)	Amount of blood present in the left ventricle at the end of diastole	50-90 mL/m <sup>2</sup>
End-systolic volume (ESV)	Amount of blood present in the left ventricle at the end of systole	25 mL/m <sup>2</sup>
Stroke volume (SV)	Amount of blood ejected from the heart in one contraction (SV = EDV - ESV)	45 ± 12 mL/m <sup>2</sup>

**TABLE 7** Pressures and volumes used in cardiac monitoring—cont'd**Volumes**

Ejection fraction (EF)	Proportion (fraction) of EDV ejected from the left ventricle during systole (EF = SV/EDV)	0.67 ± 0.07
Cardiac output (CO)	Amount of blood ejected by the heart in 1 minute	3-6 L/min
Cardiac index (CI)	Amount of blood ejected by the heart in 1 minute per square meter of body surface area (CI = CO/body surface area)	2.8-4.2 L/min/m <sup>2</sup> for a patient with 1.5 m <sup>2</sup> of body surface area

- To perform dilation of stenotic coronary arteries (angioplasty), place coronary artery stents, or perform laser atherectomy


In right-sided heart catheterization, usually the jugular, subclavian, brachial, or femoral vein is used for vascular access. In left-sided heart catheterization, usually the right femoral artery is cannulated; alternatively, however, the radial or brachial artery may be chosen (Figure 8). As the catheter is placed into the great vessels of the heart chamber, pressures are monitored and recorded. Blood samples for analysis of O<sub>2</sub> content are also obtained. After pressures are obtained, angiographic visualization of the heart chambers, valves, and coronary arteries is achieved with the injection of radiographic dye.

*Percutaneous transluminal coronary angioplasty* and *intra-coronary stents* are therapeutic procedures that can be performed during cardiac catheterization in some medical facilities.

Likewise, *atherectomy* of coronary arterial plaques can be performed to more permanently open some of the hard, atheromatous plaques.

**Contraindications**

- Patients with an iodine dye allergy who have not received preventive medication for allergy

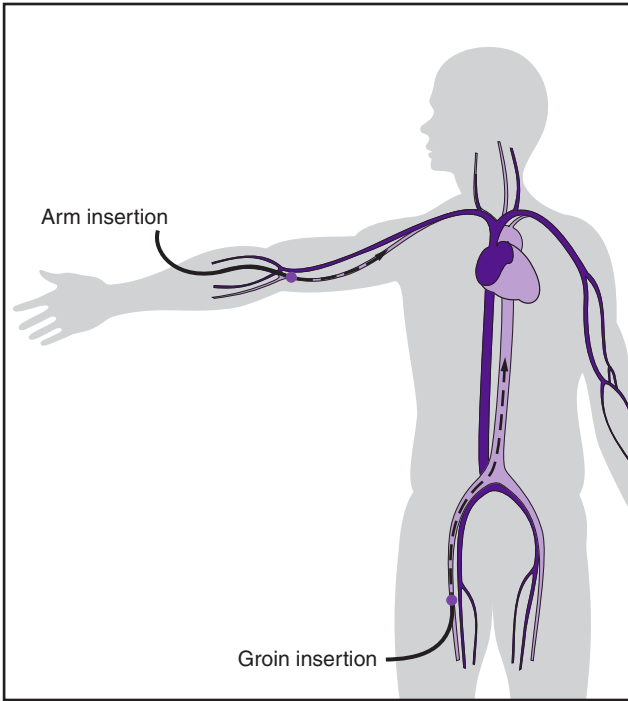

**FIG. 8** Cardiac catheterization. Brachial (arm) or femoral (groin) arterial insertion for cardiac catheterization.

- Patients who are pregnant because there is a risk of radiation exposure to fetuses
- Patients with renal disorders because iodinated contrast is nephrotoxic
- Patients with a bleeding propensity

### Potential complications

- Cardiac arrhythmias (dysrhythmias)
- Perforation of the heart myocardium
- Catheter-induced embolic stroke or myocardial infarction
- Complications associated with the catheter insertion site (e.g., arterial thrombosis, embolism, or pseudoaneurysm)
- Infection at the catheter insertion site
- Pneumothorax after subclavian vein catheterization of the right side of the heart



- Lactic acidosis in patients who are taking metformin and receiving iodine contrast. Metformin should be held the day of the test to prevent this complication.

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.

- Obtain written permission from the fully informed patient.

**PT** Allay the patient's fears and anxieties regarding this test. Although this test creates tremendous fear in a patient, it is performed often and complications are rare.

**PT** Instruct the patient to abstain from oral intake for at least 4 to 8 hours before the test.

- Prepare the catheter insertion site as per protocol.
- Mark the patient's peripheral pulses with a pen before catheterization. This will facilitate postcatheterization assessment.
- Provide appropriate precatheterization sedation as ordered.

**PT** Instruct the patient to void before the procedure.

- Remove all valuables and dental prostheses before transporting the patient to the catheterization laboratory.
- Obtain IV access for delivery of IV fluids and cardiac drugs if necessary.

### During

- Take the patient to the cardiac catheterization laboratory.
- Note the following procedural steps:
  1. The chosen catheter insertion site is prepared and draped in a sterile manner.
  2. The desired vessel is punctured, and the angiographic catheter is threaded through the sheath over a guidewire.
  3. When the catheter is in the desired location, the appropriate cardiac pressures and volumes are measured.
  4. Cardiac ventriculography is performed with controlled injection of contrast.
  5. Each coronary artery is catheterized. Cardiac angiography is then carried out with a controlled injection of contrast material.
  6. During the injection, x-ray images are rapidly made.
  7. The patient's vital signs must be monitored constantly during this procedure.
  8. After obtaining all the required information, the catheter is removed, and a vascular closure device may be placed.
  9. A chemical vascular closure device designed to seal the arterial puncture is often placed.

## 206 cardiac catheterization

- Note that this test is usually performed by a cardiologist in approximately 1 hour.
- PT** Tell the patient that during the injection, he or she may experience a severe hot flush. This is uncomfortable but lasts only 10 to 15 seconds.
- Note that some patients have a tendency to cough as the catheter is placed into the pulmonary artery.

### After

- Monitor the patient's vital signs to check for bleeding.
- Apply pressure to the site of vascular access.
- Keep the patient on bed rest for 4 to 8 hours to allow for complete sealing of the arterial puncture.
- Keep the affected extremity extended and immobilized with sandbags to decrease bleeding.
- Assess the puncture site for signs of bleeding, hematoma, or absence of pulse.
- Compare the patient's pulses with preprocedure baseline values.
- PT** Encourage the patient to drink fluids to maintain adequate hydration. Dehydration may be caused by the diuretic action of the dye. Monitor urinary output.
- PT** Instruct the patient to report any signs of numbness, tingling, pain, or loss of function in the involved extremity.
- See p. xxi for appropriate interventions concerning care of patients with iodine allergy.
- PT** Instruct the patient that the test will be reviewed by the cardiologist, and the results will be available in 1 or 2 days.

### Abnormal findings

Anatomic variation of the cardiac chambers and great vessels

Coronary artery occlusive disease

Coronary aneurysm

Coronary fistula

Cardiomyopathy

Ventricular aneurysm

Ventricular mural thrombi

Intracardiac tumor

Aortic root arteriosclerotic or aneurysmal disease

Anomalies in pulmonary venous return

Acquired or congenital septal defects and valvular abnormalities

Pulmonary emboli

Pulmonary hypertension

### notes

**cardiac nuclear scan** (Myocardial perfusion scan, Myocardial perfusion imaging, Myocardial scan, Cardiac scan, Heart scan, Thallium scan, MUGA scan, Isonitrite scan, Sestamibi cardiac scan, Cardiac flow studies, and Nuclear stress test)

**Type of test** Nuclear scan

### Normal findings

Heterogeneous uptake of radionuclide throughout the myocardium of the left ventricle

Left ventricular end-diastolic volume  $\leq 70$  mL

Left ventricular end-systolic volume  $\leq 25$  mL

Left ventricular ejection fraction  $> 50\%$

Right ventricular ejection fraction  $> 40\%$

Normal cardiac wall motion

No muscle wall thickening

### Test explanation and related physiology

A cardiac perfusion scan measures the coronary blood flow at rest and during exercise. It is often used to evaluate the cause of chest pain. It may be done after a coronary ischemic event to evaluate coronary patency or heart muscle function.

In this test, a radionuclide is injected intravenously into the patient. Myocardial perfusion images are then obtained while the patient is lying down under a single-photon emission computed tomography (SPECT) camera that generates a picture of the radioactivity coming from the heart. This scan can be performed at rest or with exercise such as treadmill or bicycling (*myocardial nuclear stress testing*). Medications may be administered that duplicate exercise stress testing. Vasodilators (dipyridamole, adenosine, and regadenoson) or chronotropic agents (dobutamine) are commonly used. Regadenoson is the most recent  $A_{2A}$  adenosine receptor agonist that instigates coronary vasodilatation.

The initial radioisotope used was thallium (thus the name *thallium scan*). Technetium agents such as tetrofosmin and sestamibi (isonitrite) are now more commonly used. The uptake of these agents is proportional to the myocardial coronary flow. At rest, a coronary stenosis must exceed 90% of the normal diameter before blood flow is impaired enough to see it on the perfusion scan. With exercise stress testing, however, stenosis of 50% becomes obvious. Myocardial perfusion scans can be synchronized by gating the images with the cardiac cycle and thereby allowing the visualization and evaluation of cardiac muscle function. Prior muscle injury is demonstrated by reduced muscle wall

## 208 cardiac nuclear scan

motion. Most times, nuclear myocardial scans include both perfusion and gated wall motion images.

Cardiac nuclear imaging when gated to the cardiac cycle (*multigated acquisition scan [MUGA]*, *gated blood pool scan*) can provide an accurate measure of ventricular function through the calculation of the *ventricular ejection fraction*. In this scan, the patient's red blood cells (RBCs) are tagged with technetium. Ventricular volumes can be calculated and used to accurately calculate the amount of blood that is ejected from the ventricle with each contraction (ejection fraction). Patients with cardiomyopathies (ischemic, infiltrative, inflammatory), cardiac transplant, or drug-induced cardiac muscle toxicity (from doxorubicin or Herceptin) require frequent evaluation of ventricular ejection fraction.

### Contraindications

- Patients who are uncooperative or medically unstable
- Patients with severe cardiac arrhythmia
- Patients who are pregnant (unless the benefits outweigh the risks) because there is risk of fetal exposure to radionuclide material

### Interfering factors

- Myocardial trauma
- Cardiac flow studies can be altered by excessive alterations in chest pressure (as exists with excessive crying in children).
- Recent nuclear scans (e.g., thyroid or bone scan)
- ✚ Drugs, such as long-acting nitrates, may only temporarily improve coronary perfusion and cardiac function.

### Procedure and patient care

#### Before

- PT Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT Instruct the patient that a short fasting period may be required, especially when using sestamibi or tetrofosmin.
- PT Tell the patient that the only discomfort associated with this test is the venipuncture required for injection of the radioisotope.
  - Be sure that all jewelry is removed from the chest wall.
  - Obtain a consent form if stress testing is to be performed.

#### During

- Take the patient to the nuclear medicine department. Depending on the type of nuclear myocardial scan, each scanning protocol is different.

- Note the following general procedural steps:
  1. One or more IV injection of radionuclide material is performed.
  2. Electrocardiographic (ECG) leads may be applied.
  3. Depending on the radionuclide used, scanning is performed 15 minutes to 4 hours later.
  4. The SPECT camera is placed at the level of the precordium.
  5. If a single gamma camera is used, the patient is placed in a supine position and then may be repositioned to the lateral position or in the right and left oblique positions. In some departments, the detector can be rotated around the patient, who remains in the supine position.
  6. The gamma ray scanner records the image of the heart, and an image is immediately developed.
  7. For an exercise stress test, additional radionuclide is injected during exercise when the patient reaches a maximum heart rate. The patient then lies on a table, and scanning is done. A repeat scan may be done 3 to 4 hours later.
- Note that myocardial scans are usually performed in less than 30 minutes by a nuclear medicine technician.
- If nuclear cardiac stress testing is performed, follow the routine protocol described on p. 212.

### After

- PT** Inform the patient that because only tracer doses of radioisotopes are used, no precautions need to be taken against radioactive exposure to personnel or family.
- PT** Instruct the patient to drink fluids to aid in the excretion of the radioactive substance.
  - Apply pressure or a pressure dressing to the venipuncture site.
  - Assess the venipuncture site for bleeding.
  - If stress testing was performed, evaluate the patient's vital signs at frequent intervals (as indicated).
  - Remove any applied ECG leads.

### Abnormal findings

Coronary artery occlusive disease

Decreased myocardial function associated with ischemia, myocarditis, cardiomyopathy, or congestive heart failure

Decreased cardiac output

### notes

**cardiac stress testing** (Exercise stress testing; Nuclear stress testing; Echo stress testing)

**Type of test** Electrodiagnostic; nuclear

**Normal findings**

Patient able to obtain and maintain maximal heart rate of 85% for predicted age and gender with no cardiac symptoms or ECG change. No cardiac muscle wall dysfunction present.

**Test explanation and related physiology**

Stress testing is used in the following situations:

- To evaluate chest pain in a patient suspected of having coronary disease
- To determine the limits of safe exercise during a cardiac rehabilitation program or to assist patients with cardiac disease in maintaining good physical fitness
- To detect labile or exercise-related hypertension
- To detect intermittent claudication in patients with suspected vascular occlusive disease in the extremities
- To evaluate the effectiveness of treatment in patients who take antianginal or antiarrhythmic medications
- To evaluate the effectiveness of cardiac intervention (e.g., bypass grafting or angioplasty)

Stress testing is a noninvasive study that provides information about the patient's cardiac function. In stress testing, the heart is stressed in some way and then evaluated during the stress. Changes indicating ischemia suggest coronary occlusive disease. By far the most commonly used method is *exercise stress testing*. *Chemical stress testing* methods are becoming more commonly used because of their safety and increased accuracy. A third method, less commonly used, is *pacemaker stress testing*.

During *exercise stress testing*, the ECG, heart rate, and blood pressure are monitored while the patient engages in some type of physical activity (stress). The treadmill test is the most frequently used because it is the most easily standardized and reproducible.

The usual goal of exercise stress testing is to increase the heart rate to just below maximal levels or to the *target heart rate*. Usually this target heart rate is 80% to 90% of the maximal heart rate. The test is usually discontinued if the patient reaches that target heart rate or develops any symptoms or ECG changes. The maximal heart rate is determined by a chart that takes into account the patient's age and gender. (Target rate is about 220 minus the patient's age.) Patients taking calcium channel

blockers and sympathetic blockers have a lower-than-expected maximal heart rate.

Exercise stress testing is based on the principle that occluded arteries will be unable to meet the heart's increased demand for blood during the testing. This may become obvious with symptoms (e.g., chest pain, fatigue, dyspnea, tachycardia, cardiac arrhythmias [dysrhythmias], fall in blood pressure) or ECG changes (e.g., ST-segment variance  $> 1$  mm, increasing premature ventricular contractions, or other rhythm disturbances). Besides the electrodiagnostic method of cardiac evaluation, the stressed heart can also be evaluated by nuclear scanning or echocardiography (which are more sensitive and accurate).

When exercise testing is not advisable or the patient is unable to exercise at a level adequate to stress the heart (e.g., patients with an orthopedic, arthritic, neurologic, vascular, or pulmonary limitation), *chemical stress testing* is recommended. Although chemical stress testing is less physiologic than exercise testing, it is safer and more controllable. *Dipyridamole* (*Persantine*) is a coronary vasodilator. If one coronary artery is significantly occluded, the coronary blood flow is diverted to the opened vessels. *Adenosine* works similarly to dipyridamole. *Dobutamine* is another chemical that can stress the heart. Dobutamine stimulates the heart muscle function. The normal heart muscle increases (augments) its contractility (wall motion). Ischemic muscle has no augmentation. In fact, in time the ischemic area becomes hypokinetic. Infarcted tissue is akinetic. In chemical stress testing, the stressed heart is evaluated by nuclear scanning or echocardiography.

*Pacing* is another method of stress testing. In patients with permanent pacemakers, the rate of capture can be increased to a rate that would be considered a cardiac stress. The heart is then evaluated electrodiagnostically or with nuclear scanning or echocardiography.

### Contraindications

- Patients with unstable angina
- Patients with severe aortic valvular heart disease
- Patients who have recently had a myocardial infarction (however, limited stress testing can be done)
- Patients with severe congestive heart failure
- Patients with severe left main coronary artery disease

### Potential complications

- Fatal cardiac arrhythmias
- Severe angina

## 212 cardiac stress testing

- Myocardial infarction
- Fainting

### Interfering factors

- Heavy meals before testing can divert blood to the GI tract.
  - Nicotine from smoking can cause coronary artery spasm.
  - Caffeine blocks the effect of dipyridamole.
  - Medical problems, such as left ventricular hypertrophy, hypertension, valvular heart disease, left bundle-branch block, severe anemia, hypoxemia, and chronic pulmonary disease, can affect results.
- ✚ Drugs that can affect test results include beta-blockers, calcium channel blockers, digoxin, and nitroglycerin.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to abstain from eating, drinking, and smoking for 4 hours before testing.
- PT** Inform the patient about the risks of the test and obtain informed consent.
- PT** Instruct the patient to bring comfortable clothing and athletic shoes for exercise. Slippers are not acceptable.
- PT** Inform the patient if any medications should be discontinued for a time period before testing.
  - Obtain a pretest ECG.
  - Record the patient's vital signs for baseline values.
  - Apply and secure appropriate ECG electrodes.

#### During

- Note that a physician usually is present during stress testing.
- After the patient begins to exercise, adjust the treadmill machine settings to apply increasing levels of stress.
- PT** Encourage patients to verbalize any symptoms.
  - Note that during the test the ECG tracing and vital signs are monitored continuously.
  - Terminate the test if the patient complains of chest pain, exhaustion, dyspnea, fatigue, or dizziness.
  - Note that testing usually takes approximately 45 minutes.
- PT** Inform the patient that the physician usually interprets and explains the results.

#### After

- Place the patient in the supine position to rest after the test.



- Monitor the ECG tracing and record vital signs at poststress intervals until recordings and values return to pretest levels.
- Remove electrodes and paste.

**PT** Tell the patient when the test results will be available.

### **Abnormal findings**

Coronary artery occlusive disease

Exercise-related hypertension or hypotension

Intermittent claudication

Abnormal cardiac rhythms: stress induced

Arrhythmias (e.g., ventricular tachycardia or supraventricular tachycardia)

**notes**

**carotid artery duplex scanning** (Carotid ultrasound)

**Type of test** Ultrasound

**Normal findings**

Carotid artery free of plaques and stenosis

**Test explanation and related physiology**

Carotid duplex scanning is a noninvasive ultrasound test used on the vertebral and extracranial carotid artery to detect occlusive disease directly. It is recommended for patients with peripheral vascular disease, headaches, and neurologic symptoms (e.g., transient ischemic attacks [TIAs], hemiparesis, paresthesia, and acute speech or visual deficits).

This scan is called *duplex* because it combines the benefits of two methods of ultrasonography—Doppler and B-mode. With the use of the transducer, a B-mode ultrasound grayscale image of the carotid vessel is obtained. A pulsed Doppler probe within the transducer is used to evaluate blood flow velocity and direction in the artery and to measure the amplitude and waveform of the carotid arterial pulse. A computer combines that information and provides a two-dimensional image of the carotid artery along with an image of blood flow. With this technique, one is able to directly visualize areas of stenotic or occluded arteries and arterial flow disruption. The degree of occlusion is measured in the percentage of the entire lumen that is occluded.

*Color Doppler ultrasound (CDU)* can be added to duplex scanning. CDU assigns color for direction of blood flow within the vessel, and the intensity of that color is dependent on the mean computed velocity of blood traveling in the vessel. This allows visualization of stenotic areas by seeing slowing or reversal of direction of blood flow at a particular area of the artery.

Measurement of the thickness of the wall of the carotid artery (*carotid intima-media thickness [CIMT]*) is used as a measurement of cerebrovascular atherosclerosis specifically and is a predictor of coronary atherosclerosis in general. CIMT is also used to monitor progression of atherosclerosis (particularly in diabetics). It is used to monitor atherosclerotic regression in patients who are undergoing a treatment for atherosclerosis.

Transcranial Doppler ultrasonography is useful for evaluating more proximal vascular anatomy, including the middle cerebral artery (MCA), intracranial carotid artery, and vertebrobasilar artery.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no special preparation is required.
- PT** Assure the patient that the test is painless.

### During

- Place the patient in the supine position with the head supported to prevent lateral motion.
- Note the following procedural steps:
  1. A water-soluble gel is used to couple the sound from the transducer to the skin surface.
  2. Images of the carotid artery and pulse waveform are obtained.
- Note that this test is performed by an ultrasound technologist in the ultrasound or radiology department in approximately 15 to 30 minutes.
- PT** Tell the patient that no discomfort is associated with this test.

### After

- Remove the water-soluble gel from the patient.

## Abnormal findings

Carotid artery occlusive disease

Carotid artery aneurysm

notes

**cell culture drug resistance testing** (CCDRT, Chemosensitivity assay, Drug response assay)

**Type of test** Miscellaneous

**Normal findings**

Cells sensitive to planned therapeutic drugs

**Test explanation and related physiology**

Cell culture drug resistance testing (CCDRT) refers to testing the reaction of a patient's own cancer cells in the laboratory to drugs that may be used to treat the patient's cancer. The idea is to identify which drugs are more likely to work and which drugs are less likely to work. By avoiding the latter and choosing from among the former, the patient's probability of benefiting from the chemotherapy may be improved.

Multiple tests are available for drug sensitivity testing, but all have four common steps. Cancer cells from the patient's tumor must be obtained and isolated. The cells are then incubated with various potentially therapeutic drugs. Assessment of cell survival is then performed, and the results are provided. Based on those results, the clinician can recommend more appropriate chemotherapy for a particular cancer. In most cases, this testing is used for patients with refractory or recurrent epithelial tumors (usually breast or ovarian cancer).

**Procedure and patient care****Before**

**pr** Explain the process to the patient. (Tumor cells are usually obtained by a surgical procedure.)

**During**

- Tumor cells are sent to a reference laboratory. The method of tissue preservation varies among laboratories.

**After**

- After the results are obtained, appropriate chemotherapy targeted to the patient's tumor cells is administered.

**Abnormal findings**

- Epithelial cancer

notes

**cell-free maternal DNA testing** (Prenatal cell-free DNA screening, Noninvasive prenatal testing [NIPT], cell-free DNA in maternal blood, cf DNA screening)

**Type of test** Blood

### Normal findings

Low risk of chromosomal abnormality

### Test explanation and related physiology

Trisomy 21 (Down syndrome), trisomy 18 (Edwards syndrome), and trisomy 13 (Patau syndrome) are the three most common chromosomal abnormalities affecting live births. Although 1 in 450 live births has one of these aneuploidy abnormalities, trisomy 21 is the most common. Abnormal findings on pelvic ultrasonography of the fetus, including fetal nuchal translucency or thickness (see p. 685) along with biochemical markers (e.g., hCG, see p. 515; PAAP-A, see p. 730), can identify pregnancies at high risk for these chromosomal defects. The definitive diagnosis requires chorionic villus sampling (CVS; see p. 241) and amniocentesis (see p. 49), which are invasive and increase the risk for miscarriage.

Cell-free (cf) DNA from the placental fetal cells circulates in maternal blood. This DNA can be extracted and, through advanced laboratory techniques of targeting genomic sequencing, allows 99% of the cases of trisomy to be detected. False-positive rates instigating unnecessary invasive testing are less than 1%. This testing can be performed as early as 10 weeks of gestation but is typically done between 10 and 22 weeks. Because of newer laboratory techniques of multiplexing, results can be available in about 1 week.

In 2012, the American College of Obstetricians and Gynecologists (ACOG) and the Society for Maternal-Fetal Medicine (SMFM) issued a joint committee opinion that supported noninvasive prenatal testing that uses cell-free fetal DNA for women at increased risk for having a baby with a chromosomal abnormality (Box 4).

Testing of high-risk pregnant women may be done in several ways:

- cf DNA test done in the first trimester along with ultrasonography
- cf DNA test done in the second trimester without ultrasonography
- A combination of both

**BOX 4 Women at high risk for having children with chromosomal abnormalities**

- Maternal age 35 years or older at delivery
- Fetal ultrasonographic findings indicating an increased risk for trisomy
- History of pregnancy with trisomy
- Positive maternal screen result
- Other translocation abnormalities with increased risk of trisomy

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
  - Fasting: no
  - Blood tube commonly used: red
- PT** Encourage all women undergoing cf DNA to have genetic counseling.
- PT** Results should be reviewed with the patient, and the risks, benefits, and alternatives to further testing should be explained.

**Abnormal findings**

Trisomy 21 (Down syndrome)  
Trisomy 18 (Edwards syndrome)  
Trisomy 13 (Patau syndrome)  
Duchenne muscular dystrophy  
Hemophilia

**notes**

**cell surface immunophenotyping** (Flow cytometry cell surface immunophenotyping, Lymphocyte immunophenotyping, AIDS T-lymphocyte cell markers, CD4 marker, CD4/CD8 ratio, CD4 percentage)

**Type of test** Blood

### Normal findings

Cells	Percent	No. of cells/ $\mu$ L
T-cells	60-95	800-2500
T-helper (CD4) cells	60-75	600-1500
CD8 T-cells	25-30	300-1000
B-cells	4-25	100-450
Natural killer cells	4-30	75-500
<i>CD4/CD8 ratio</i> : > 1		

### Test explanation and related physiology

This test is used to detect the progressive depletion of CD4 T-lymphocytes, which is associated with an increased likelihood of clinical complications from acquired immunodeficiency syndrome (AIDS). Test results can indicate if a patient with AIDS is at risk for developing opportunistic infections. It is also used to confirm the diagnosis of acute myelocytic leukemia (AML) and to differentiate AML from acute lymphocytic leukemia (ALL).

CD4 (T-helper cells) and CD8 (T-suppressor cells) are examples of T-lymphocytes. T-lymphocytes, and especially CD4 counts, when combined with HIV viral load testing (see p. 503) are, used to determine the time to initiate antiviral therapy. They also can be used to monitor antiviral therapy. Successful antiviral therapy is associated with an increase in CD4 counts. Worsening of disease or unsuccessful therapy is associated with decreasing T-lymphocyte counts.

There are three related measurements of CD4 T-lymphocytes. The first measurement is the *total CD4 count*. This is measured in whole blood and is the product of the white blood cell (WBC) count, the lymphocyte differential count, and the percentage of lymphocytes that are CD4 T-cells. The second measurement, the *CD4 percentage*, is a more accurate prognostic marker. The third prognostic marker, which is also more reliable than the total CD4 count, is the *CD4-to-CD8 ratio*.

## 220 cell surface immunophenotyping

Of the three T-cell measurements, the total CD4 count is the most variable. There is substantial diurnal variation in this count. Because it is a calculated measurement, the combination of possible laboratory error and personal fluctuation can result in wide variations in test results. With the CD4 percentage and the CD4-to-CD8 ratio, very little diurnal variation and laboratory error exist. The Multicenter AIDS Cohort Study suggests that the latter two measurements are more accurate than the total CD4 count. However, because the total CD4 count was originally thought to be the best marker, this test was used in many of the studies that now form the basis for practice recommendations.

CD4 measurement is a prognostic marker that can indicate whether a patient infected with HIV is at risk for developing opportunistic infections. The measurement of CD4 levels is used to decide whether to initiate *Pneumocystis jirovecii* pneumonia prophylaxis and antiviral therapy and for determining the prognosis of patients with HIV.

Both immunodeficiency and the dosage of immunosuppressive medications used after organ transplant are also monitored with the use of this cell surface immunophenotyping. Lymphomas and other lymphoproliferative diseases are now classified and treated according to the predominant lymphocyte type identified.

### Contraindications

- Patients who are not emotionally prepared for the prognosis that the results may indicate

### Interfering factors

- Diurnal variation occurs.
- A recent viral illness can decrease total T-lymphocyte counts.
- Nicotine and very strenuous exercise have been shown to decrease lymphocyte counts.
- ✚ Steroids can *increase* lymphocyte counts.
- ✚ Immunosuppressive drugs will *decrease* lymphocyte counts.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: green or purple
- Never recap needles. Dispose of needles and syringes in a puncture-proof container.
- Keep the specimen at room temperature. Do not refrigerate.
- The specimen must be evaluated within 24 hours.



- PT** Instruct the patient to observe the venipuncture site for infection. Patients with leukemia or AIDS are immunocompromised and susceptible to infection.
- PT** Encourage the patient to discuss his or her concerns regarding the prognostic information obtained by these results.
- Do not give test results over the phone. Decreasing CD4 counts can have devastating consequences.

### Abnormal findings

#### ▲ Increased counts

Leukemias  
Lymphoma

#### ▼ Decreased counts

Organ transplant patients  
Immunodeficiency diseases

notes

## ceramides

**Type of test** Blood

### Normal findings

Low ceramide risk score

### Test explanation and related physiology

Plasma ceramides are new biomarkers for predicting cardiovascular (CV) disease. Ceramides are complex lipids that are important in maintaining cell membrane integrity, the cellular stress response, inflammatory signaling, and cell death. Synthesis of ceramides occurs in all tissues. Elevated concentrations of ceramides are associated with atherosclerotic plaque formation, ischemic heart disease, myocardial infarction, hypertension, stroke, type 2 diabetes mellitus, and metabolic syndromes. Even after these new markers have been adjusted for age, gender, smoking status, and serum biomarkers (e.g., low-density lipoprotein [LDL] and high-density lipoprotein [HDL] cholesterol, C-reactive protein [CRP], and Lp-PLA2), ceramides continue to be associated with the risk of CV disease. Therefore plasma ceramides are predictors of future adverse CV events as described in the table below in which a risk score is based on the absolute levels of the ceramides and analytic mathematical divisions of those levels.

#### Ceramides as Predictors of Adverse CV Events

Ceramide score	Relative risk	Risk category
0-2	1.0	Lower
3-6	1.8	Moderate
7-9	2.3	Increased
10-12	5.1	Higher

Three specific ceramides have been identified as highly linked to CV disease and insulin resistance: Cer16:0, Cer18:0, and Cer24:1. Plasma ceramides can predict adverse CV events within 1 to 5 years depending on existence of other cardiac disease. Plasma ceramides can be reduced by treatment with simvastatin, rosuvastatin, and ezetimibe. Because these scores improve within a few weeks of treatment, knowing their ceramide scores may improve patient motivation and therapeutic compliance.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

**Abnormal findings**

Increased risk for CV event.

notes

**cervical biopsy** (LEEP procedure, Cone biopsy)

**Type of test** Microscopic examination

**Normal findings**

Normal squamous cells

**Possible critical values**

Cancer cells

**Test explanation and related physiology**

When a Pap smear reveals an *epithelial cell abnormality* or when a pelvic examination reveals a possible abnormality in the cervix, a biopsy of that structure is performed. There are several different methods of biopsy, all of which obtain an increasing amount of tissue. Cervical biopsy procedures include the following:

A *simple cervical biopsy*, sometimes called a *punch biopsy*, removes a small piece of tissue from the surface of the cervix. This is often performed during colposcopy (see p. 262).

An *endocervical biopsy* (*endocervical curettage*) removes tissue from high in the cervical canal by scraping with a sharp instrument.

*Loop electrosurgical excision procedure* (LEEP) uses a thin, low-voltage electrified wire loop to cut out abnormal tissue on the cervix and high in the endocervical canal (sometimes called a *large loop excision of the transformation zone* [LLETZ]).

A *cone biopsy* (*conization*) is a more extensive form of a cervical biopsy. It is called a cone biopsy because a cone-shaped wedge of tissue is removed from the cervix. Both normal and abnormal cervical tissues are removed. This can be performed by LEEP, surgical knife (scalpel), or carbon dioxide laser.

**Contraindications**

- Patients with active menstrual bleeding
- Patients who are pregnant

**Potential complications**

- After the surgery, a small number of women (< 10%) may have significant bleeding that requires vaginal packing or a blood transfusion.
- Infection of the cervix or uterus may occur. (This is rare.)
- Narrowing of the cervix (cervical stenosis) may occur, which can cause infertility. (This is rare.)

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- Obtain informed consent if required by the institution.

### During

- Note the following procedural steps:
    1. The patient is placed in the lithotomy position, and a vaginal speculum is used to expose the vagina and cervix.
    2. The cervix is cleansed with a 3% acetic acid solution or antiseptic to remove excess mucus and cellular debris and to accentuate the difference between normal and abnormal epithelial tissues.
    3. Medication is injected to numb the cervix (*cervical block*).
    4. With the instrument chosen by the doctor, a punch biopsy, endocervical biopsy, LEEP, or cone biopsy is performed.
  - Note that the physician performs the procedure in approximately 30 minutes.
  - Whereas cone biopsy is done in the operating room, the other procedures can be performed in the doctor's office.
- PT** Tell the patient that some women complain of pressure pains from the vaginal speculum and that discomfort may be felt if biopsy specimens are obtained.
- Most women can return to normal activities immediately after a simple cervical biopsy or an endocervical biopsy.
  - Most women will be able to return to normal activities within 2 to 4 days after LEEP or cone biopsies. This can vary, depending on the amount of tissue removed.

### After

- PT** Inform the patient that it is normal to experience the following:
1. Vaginal bleeding if biopsy specimens were taken. Suggest that she wear a sanitary pad.
  2. Mild cramping for several hours after the procedure
  3. Brownish-black vaginal discharge during the first week
  4. Vaginal discharge or spotting for about 1 to 3 weeks
- PT** Instruct the patient to use sanitary napkins instead of tampons for 1 to 3 weeks.
- PT** Tell the patient to avoid sexual intercourse for 3 to 4 weeks.
- PT** Inform the patient not to douche for 3 to 4 weeks.
- PT** Tell the patient how to obtain the test results.

## 226 cervical biopsy

**PT** Instruct the patient to call the doctor for any of the following symptoms:

1. Fever
2. Spotting or bleeding that lasts longer than 1 week
3. Bleeding that is heavier than a normal menstrual period and contains blood clots
4. Increasing pelvic pain
5. Bad-smelling, yellowish vaginal discharge, which may indicate an infection

### **Abnormal findings**

Chronic cervical infection

Cervical intraepithelial neoplasia

Cervical carcinoma *in situ*

Invasive cervical carcinoma

Endocervical adenocarcinoma

### **notes**

## chest x-ray (CXR)

**Type of test** X-ray

### Normal findings

Normal lungs and surrounding structures

### Test explanation and related physiology

The CXR image is important in a complete evaluation of the pulmonary and cardiac systems. Much information can be provided by CXR images. One can identify or follow (by repeated CXR images) the following:

- Tumors of the lung (primary and metastatic), heart (myxoma), chest wall (soft tissue sarcomas), and bony thorax (osteogenic sarcoma)
- Inflammation of the lung (pneumonia), pleura (pleuritis), and pericardium (pericarditis)
- Fluid accumulation in the pleura (pleural effusion), pericardium (pericardial effusion), and lung (pulmonary edema)
- Air accumulation in the lung (chronic obstructive pulmonary disease) and pleura (pneumothorax)
- Fractures of the bones of the thorax or vertebrae
- Diaphragmatic hernia
- Heart size, which may vary depending on cardiac function
- Calcification, which may indicate large-vessel deterioration or old lung granulomas
- Location of centrally placed IV access devices

Most CXRs are taken with the patient standing. The sitting or supine position also can be used, but x-ray images taken with the patient in the supine position will not demonstrate fluid levels. A *posteroanterior (PA)* view, with the x-rays passing through the back of the body (posterior) to the front of the body (anterior), is taken first. Then a *lateral* view, with the x-rays passing through the patient's side, is taken.

*Oblique* views may be taken with the patient turned at different angles as the x-rays pass through the body. *Lordotic* views provide visualization of the apices (rounded upper portions) of the lungs and are usually used for detection of tuberculosis. *Decubitus* images are taken with the patient in the recumbent lateral position to localize fluid, which becomes dependent within the pleural space (pleural effusion).

Chest x-ray studies are best performed in the radiology department. Studies using a portable x-ray machine may be done

## 228 chest x-ray

at the bedside and are often performed on critically ill patients who cannot leave the nursing unit.

### Contraindications

- Patients who are pregnant unless the benefits outweigh the risks

### Interfering factors

- Conditions (e.g., severe pain) that prevent the patient from taking and holding a deep breath
- Scarring from previous lung surgery, which makes interpretation difficult
- Obesity, which requires more x-rays to penetrate the body to provide a readable picture

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Tell the patient that no fasting is required.
- PT** Instruct the patient to remove clothing to the waist and to put on an x-ray gown.
- PT** Inform the patient to remove all metal objects (e.g., necklaces, pins) so that they do not block visualization of part of the chest.
- PT** Tell the patient that he or she will be asked to take a deep breath and hold it while the x-ray images are taken.
- PT** Instruct men to ensure that their testicles are covered and women to have their ovaries covered, using a lead shield to prevent radiation-induced abnormalities.

#### During

- PT** After the patient is correctly positioned, tell him or her to take a deep breath and hold it until the images are taken.
- Note that x-ray images are taken by a radiologic technologist in several minutes.
- PT** Inform the patient that no discomfort is associated with chest radiography.

#### After

- Note that no special care is required after the procedure.



**Abnormal findings****Lung**

Lung tumor (primary or metastatic)  
 Pneumonia  
 Pulmonary edema  
 Pleural effusion  
 Chronic obstructive pulmonary disease  
 Pneumothorax  
 Atelectasis  
 Tuberculosis  
 Lung abscess  
 Congenital lung diseases (hypoplasia)  
 Pleuritis  
 Foreign bodies (chest, bronchus, or esophagus)

**Heart**

Cardiac enlargement  
 Pericarditis  
 Pericardial effusion

**notes****Chest wall**

Soft tissue sarcoma  
 Osteogenic sarcoma  
 Fracture (ribs or thoracic spine)  
 Thoracic spine scoliosis  
 Metastatic tumor to the bony thorax

**Diaphragm**

Diaphragmatic or hiatal hernia

**Mediastinum**

Aortic calcinosis  
 Enlarged lymph nodes  
 Dilated aorta  
 Thymoma  
 Lymphoma  
 Substernal thyroid  
 Widened mediastinum

## Chlamydia

**Type of test** Blood; microscopic examination

### Normal findings

Negative culture

Antibodies:

*Chlamydophila pneumoniae*

IgG: < 1:64

IgM: < 1:10

*Chlamydophila psittaci*

IgG: < 1:64

IgM: < 1:10

*Chlamydia trachomatis*

IgG: < 1:64

IgM: < 1:10

Nucleic acid detection: negative

### Test explanation and related physiology

There are many *Chlamydia* species that cause various diseases within the human body. *Chlamydophila psittaci* causes respiratory tract infections, headache, altered mentation, and hepatosplenomegaly. It occurs as a result of close contact with infected birds. *Chlamydophila pneumoniae*, another species, causes pneumonia. *Chlamydophila trachomatis* infection is probably the most frequently occurring sexually transmitted disease in developed countries. Infections of the genitalia, pelvic inflammatory disease, urethritis, cervicitis, salpingitis, and endometritis are most common. *C. trachomatis* may also infect the conjunctiva, pharynx, urethra, and rectum and cause lymphogranuloma venereum. The second serotype of *C. trachomatis* causes the eye disease *trachoma*, which is the most common form of preventable blindness. A third serotype produces genital and urethral infections different from lymphogranuloma. Most women colonized with *Chlamydia* are asymptomatic.

Tests can be performed on the blood of infected patients or swabs from the conjunctiva, nasopharynx, urethra, rectum, vagina, or cervix. Urine, seminal fluid, or pelvic washing can be used in culture and in direct identification of *Chlamydia*. Early identification of infection enables sexual partners to seek testing and/or treatment as soon as possible and reduces the risk of disease spread. Prompt treatment reduces the risk of infertility in women.

## Interfering factors

- Women presently having routine menses
- Patients undergoing antibiotic therapy

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient.

### During

- For *Chlamydia* antibody testing, collect venous blood in a red-top tube. Acute and convalescent serum should be drawn 2 to 3 weeks apart.
- Sputum cultures (see p. 841) are used to check for *C. psittaci* respiratory infections.
- *Chlamydia* can be tested by direct nucleic acid identification or culture. A conjunctival smear is obtained by swabbing the eye lesion with a cotton-tipped applicator or scraping with a sterile ophthalmic spatula and smearing on a clean glass slide.
- Note the following procedural steps for *cervical culture*:
  1. The patient should refrain from douching and bathing in a tub before the cervical culture is performed.
  2. The patient is placed in the lithotomy position.
  3. A nonlubricated vaginal speculum is inserted to expose the cervix.
  4. Excess mucus is removed using a cleaning swab.
  5. A second sterile, cotton-tipped swab is inserted into the endocervical canal and moved from side to side for 30 seconds to obtain the culture.
  6. The swab is then placed into an appropriate transport tube for testing.
- Note the following procedural steps for *urethral culture*:
  1. The urethral specimen should be obtained from the man before voiding within the previous hour.
  2. A culture is taken by inserting a thin sterile swab with rotating movement about 3 to 4 cm into the urethra.
- Note the following procedural steps for a *urine specimen*:
  1. The patient should not have urinated for at least 1 hour before specimen collection.
  2. The patient should collect the first portion (first part of stream) of a random voided urine into a sterile, plastic, preservative-free container.
  3. Transfer 2 mL of urine into the urine specimen collection tube using the disposable pipette provided. (The correct volume of urine has been added when the fluid level is between the black fill lines on the urine tube.)

## 232 *Chlamydia*

- Note that these tests are performed by a physician, nurse, or other health-care provider in several minutes.
- PT** Tell the patient that these procedures cause minimal discomfort.

### **After**

- Treat patients who have positive smears with antibiotics.
- PT** Tell affected patients to have their sexual partners examined.

### **Abnormal findings**

*Chlamydia* infection

notes

## chloride, blood (Cl)

**Type of test** Blood

### Normal findings

Adult/elderly: 98-106 mEq/L or 98-106 mmol/L (SI units)

Child: 90-110 mEq/L

Newborn: 96-106 mEq/L

Premature infant: 95-110 mEq/L

### Possible critical values

< 80 or > 115 mEq/L

### Test explanation and related physiology

This test is performed as a part of multiphasic testing in what is usually called electrolytes, Astra-7, chemprofile, metabolic panel, kidney profile, comprehensive profile, SMA 12, or SMA 6. By itself, not much information is obtained. However, with interpretation of the other electrolytes, chloride can give an indication of acid-base balance and hydrational status.

Hypochloremia and hyperchloremia rarely occur alone and are usually part of parallel shifts in sodium or bicarbonate levels. Signs and symptoms of hypochloremia include hyperexcitability of the nervous system and muscles, shallow breathing, hypotension, and tetany. Signs and symptoms of hyperchloremia include lethargy, weakness, and deep breathing.

### Interfering factors

- Excessive infusions of saline can result in increased chloride levels.
- ✚ Drugs that may cause *increased* serum chloride levels include acetazolamide, ammonium chloride, androgens, chlorothiazide, cortisone preparations, estrogens, guanethidine, hydrochlorothiazide, methyldopa, and nonsteroidal antiinflammatory drugs (NSAIDs).
- ✚ Drugs that may cause *decreased* levels include aldosterone, bicarbonates, corticosteroids, cortisone, hydrocortisone, loop diuretics, thiazide diuretics, and triamterene.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or green

### Abnormal findings

#### ▲ Increased levels (hyperchloremia)

- Dehydration
- Renal tubular acidosis
- Excessive infusion of normal saline
- Cushing syndrome
- Eclampsia
- Multiple myeloma
- Kidney dysfunction
- Metabolic acidosis
- Hyperventilation
- Anemia
- Respiratory alkalosis
- Hyperparathyroidism

#### ▼ Decreased levels (hypochloremia)

- Overhydration
- Congestive heart failure
- Syndrome of inappropriate antidiuretic hormone
- Vomiting
- Chronic gastric suction
- Chronic respiratory acidosis
- Salt-losing nephritis
- Addison disease
- Burns
- Metabolic alkalosis
- Diuretic therapy
- Hypokalemia
- Aldosteronism
- Respiratory acidosis

notes

## cholesterol

**Type of test** Blood

### Normal findings

Adult/elderly: < 200 mg/dL or < 5.2 mmol/L (SI units)

Child: 120-200 mg/dL

Infant: 70-175 mg/dL

Newborn: 53-135 mg/dL

### Test explanation and related physiology

Cholesterol is the main lipid associated with arteriosclerotic vascular disease. Nearly 75% of the cholesterol is bound to LDLs and 25% is bound to HDLs. Therefore cholesterol is the main component of LDLs and only a minimal component of HDLs and very low-density lipoproteins. LDLs are most directly associated with increased risk of coronary heart disease (CHD).

The purpose of cholesterol testing is to identify patients at risk for arteriosclerotic heart disease. Cholesterol testing is usually done as a part of *lipid profile* testing, which also evaluates lipoproteins (see p. 565) and triglycerides (see p. 908), because by itself cholesterol is not a totally accurate predictor of heart disease.

Because the liver is required to make cholesterol, low serum cholesterol levels are indicative of severe liver diseases. Furthermore, because our main source of cholesterol is our diet, malnutrition is also associated with low cholesterol levels. Certain illnesses can affect cholesterol levels. For example, patients with an acute myocardial infarction may have as much as a 50% reduction in cholesterol level for as many as 6 to 8 weeks.

The goal for high-risk patients (those with known coronary artery disease or more than two risk factors) is LDL less than 70 mg/dL. Many prospective studies have shown that high serum concentrations of LDL-C are a major risk factor for CHD. Moreover, lowering of LDL-C levels will reduce the risk for major coronary events. The cholesterol-to-HDL ratio has been used to assess the risk of CHD (Table 8). Familial hyperlipidemias and hyperlipoproteinemias are often associated with high cholesterol.

### Interfering factors

- Pregnancy is usually associated with elevated cholesterol levels.
- Oophorectomy increases levels.

**TABLE 8** Cholesterol-to-HDL ratio as an indicator of risk of coronary heart disease

Risk	Ratio	
	Male	Female
½ Average	3.4	3.3
Average	5	4.4
2 × Average	10	7
3 × Average	24	11

*HDL*, High-density lipoprotein.

- ▼ Drugs that may cause *increased* levels include adrenocorticotropic hormone, anabolic steroids, beta-adrenergic blocking agents, corticosteroids, cyclosporine, epinephrine, oral contraceptives, phenytoin (Dilantin), sulfonamides, thiazide diuretics, and vitamin D.
- ▼ Drugs that may cause *decreased* levels include allopurinol, androgens, bile salt-binding agents, captopril, chlorpropamide, clofibrate, colchicine, colestipol, erythromycin, isoniazid, liothyronine (Cytomel), monoamine oxidase inhibitors, neomycin (oral), niacin, nitrates, and statins.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red
- PT** Instruct the patient to fast 12 to 14 hours after eating a low-fat diet before testing. Only water is permitted. Food can elevate triglyceride levels.
- PT** Inform the patient that dietary intake at least 2 weeks before testing will affect results.
- PT** Tell the patient that no alcohol should be consumed 24 hours before the test.
- The fingerstick method is often used in mass screening.
- PT** Instruct patients with high levels regarding a low-cholesterol diet, exercise, and appropriate body weight.



## Abnormal findings

### ▲ Increased levels

- Hypercholesterolemia
- Hyperlipidemia
- Hypothyroidism
- Uncontrolled diabetes mellitus
- Nephrotic syndrome
- Pregnancy
- High-cholesterol diet
- Xanthomatosis
- Hypertension
- Myocardial infarction
- Atherosclerosis
- Biliary cirrhosis
- Stress
- Nephrosis

### ▼ Decreased levels

- Malabsorption
- Malnutrition
- Hyperthyroidism
- Cholesterol-lowering medication
- Pernicious anemia
- Hemolytic anemia
- Sepsis
- Stress
- Liver disease
- Acute myocardial infarction

## notes

**cholinesterase** (ChE, Pseudocholinesterase [PChE], Cholinesterase RBC, Red blood cell cholinesterase, Acetylcholinesterase)

**Type of test** Blood

### Normal findings

Serum cholinesterase: 8-18 units/mL or 8-18 units/L (SI units)

RBC cholinesterase: 5-10 units/mL or 5-10 units/L (SI units)

Dibucaine inhibition: 79%-84%

(Values vary with laboratory test methods.)

### Test explanation and related physiology

This test is done to identify patients with pseudocholinesterase deficiency before anesthesia or to identify patients who may have been exposed to phosphate poisoning. Cholinesterases hydrolyze acetylcholine and other choline esters and thereby regulate nerve impulse transmission at the nerve synapse and neuromuscular junction. There are two types of cholinesterases: *acetylcholinesterase*, also known as *true cholinesterase*, and pseudocholinesterase. True cholinesterase exists primarily in the RBCs and nerve tissue. It is not in the serum. *Pseudocholinesterase*, on the other hand, exists in the serum. Deficiencies in either of these enzymes can be acquired or congenital.

Because succinylcholine (the most commonly used muscle relaxant during anesthesia induction) is inactivated by pseudocholinesterase, people with an inherited pseudocholinesterase enzyme deficiency exhibit increased and/or prolonged effects of succinylcholine. Prolonged muscle paralysis and apnea will occur after anesthesia in these patients. This situation can be avoided by measuring serum cholinesterase (pseudocholinesterase) in all patients with a family history of prolonged apnea after surgery.

Because patients with a nonfunctioning variant of pseudocholinesterase will have normal total quantitative pseudocholinesterase levels yet still have prolonged paralytic effects of succinylcholine, a second test (dibucaine inhibition) usually is also performed. Dibucaine is a known local anesthetic that inhibits the function of normal pseudocholinesterase. The *dibucaine inhibition number* is the percent of pseudocholinesterase activity that is inhibited when dibucaine is added to the patient's serum sample. If total pseudocholinesterase is normal and dibucaine numbers are low, the presence of a nonfunctioning pseudocholinesterase variant is suspected, and the patient will be at risk for succinylcholine-induced prolonged paralysis.

A common form of acquired cholinesterase deficiency, either true or pseudocholinesterase, is caused by overexposure to pesticides, organophosphates, or nerve gas. The half-life of the pseudoenzyme in serum is about 8 days, and the true cholinesterase (acetylcholinesterase [AChE]) of RBCs is more than 3 months (determined by erythropoietic activity). Recent exposure up to several weeks is determined by assay of the pseudoenzyme and months after exposure by measurement of the red cell enzyme. Persons with jobs associated with chronic exposure to these chemicals are often monitored by the frequent testing of RBC cholinesterase activity. Other potential causes of reduced cholinesterase activity include chronic liver diseases, malnutrition, and hypoalbuminemia.

Increased cholinesterase activity, when found in the amniotic fluid, represents strong evidence for a *neural tube defect (NTD)*. When an NTD is suspected based upon maternal serum alpha-fetoprotein (AFP) screening results or diagnosed via ultrasound, analyses of AFP and AChE in amniotic fluid are useful diagnostic tools.

### Interfering factors

- Pregnancy decreases test values.
- ✱ Drugs that may cause *decreased* values include atropine, caffeine, codeine, estrogens, morphine sulfate, neostigmine, oral contraceptives, phenothiazines, quinidine, theophylline, and vitamin K.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- If the test is done to identify the presurgical patient who may be at risk for cholinesterase deficiency, be sure that the test is completed several days before the planned surgery.
- It may be recommended to withhold medications that could alter test results for 12 to 24 hours before the test.

## 240 cholinesterase

### Abnormal findings

#### ▲ Increased serum levels

- Hyperlipidemia
- Nephrosis
- Diabetes mellitus

#### ▲ Increased RBC levels

- Reticulocytosis
- Sickle cell disease

#### ▼ Decreased serum levels

- Poisoning from organophosphate insecticides
- Hepatocellular disease
- Persons with congenital pseudocholinesterase enzyme deficiency
- Malnutrition

#### ▼ Decreased RBC levels

- Congenital cholinesterase deficiency
- Poisoning from organophosphate insecticides

notes

## chorionic villus sampling (CVS, Chorionic villus biopsy [CVB])

**Type of test** Cell analysis

### Normal findings

No genetic or biochemical disorders

### Test explanation and related physiology

CVS is performed on women whose unborn children may be at risk for life-threatening or significant life-altering genetic defects. This includes women who:

- Are older than age 35 years at the time of pregnancy
- Have had frequent spontaneous abortions
- Have had previous pregnancies with fetuses or infants with chromosomal or genetic defects (e.g., Down syndrome)
- Have a genetic defect themselves (e.g., hemoglobinopathies)

CVS can be performed between 8 and 12 weeks of gestation for the early detection of genetic and biochemical disorders. Because CVS detects congenital defects early, first-trimester therapeutic abortions can be performed if indicated.

For this study, a sample of chorionic villi from the chorion frondosum, which is the trophoblastic origin of the placenta, is obtained for analysis. These villi are present from 8 to 12 weeks on and reflect fetal chromosome, enzyme, and deoxyribonucleic acid content, thus permitting a much earlier diagnosis of prenatal problems than amniocentesis (see p. 49), which cannot be done before 14 to 16 weeks. Furthermore, the cells derived by CVS are more easily grown in tissue culture for *karyotyping* (determination of chromosomal or genetic abnormalities). The cells obtained with amniocentesis take a longer time to grow in culture, further adding to the delay in obtaining results.

### Potential complications

- Accidental abortion
- Infection
- Bleeding
- Fetal limb deformities

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Be certain that the physician has obtained a signed consent.


**PT** Tell the patient that no food or fluid restrictions are necessary.

## 242 chorionic villus sampling

- PT** Encourage the patient to drink at least 1 to 2 glasses of fluid before testing.
- PT** Instruct the patient not to urinate for several hours before testing. A full bladder is an excellent reference point for pelvic ultrasound.
- Assess the vital signs of the mother and the fetal heart rate of the fetus before testing. These are baseline studies that should be repeated during and on completion of the test.

### During

- Note the following procedural steps:
  1. The patient is placed in the lithotomy position.
  2. Samples of vaginal mucus may be obtained to rule out pre-procedural infections (e.g., *Chlamydia*).
  3. A cannula from the endoscope is inserted into the cervix and uterine cavity (Figure 9).
  4. Under ultrasound guidance, the cannula is rotated to the site of the developing placenta.

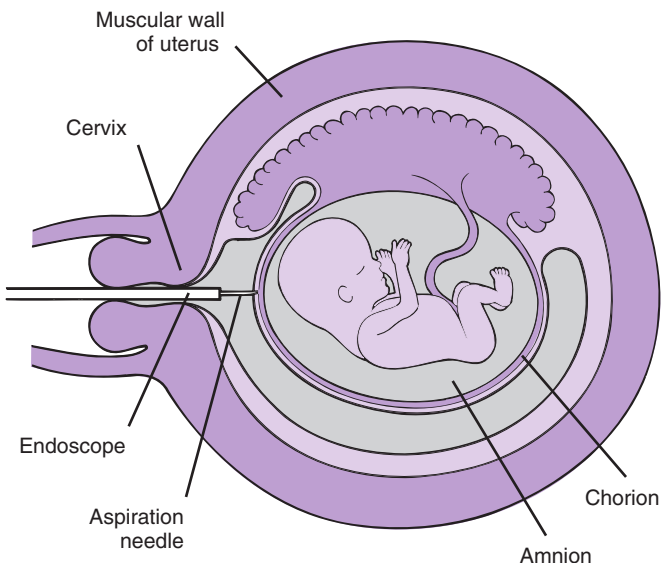

**FIG. 9** Chorionic villus sampling. Diagram of an 8-week pregnancy showing endoscopic aspiration of extraplacental villi.

5. A syringe is attached, and suction is applied to obtain several samples of villi.
  6. As many as three or more samples may be obtained to get sufficient tissue for accurate sampling.
  7. If ultrasound indicates that the trophoblastic tissue is remote from the cervix, a transabdominal approach similar to that described for amniocentesis (see p. 49) may be used.
- Note that this procedure is performed by an obstetrician in approximately 30 minutes.
- PT** Inform the patient that discomfort associated with this test is similar to that of a Pap smear.

### After

- Note that Rh-negative women (who have not been sensitized to Rh incompatibility) receive RhoGAM. This procedure may be contraindicated for women with known preexisting Rh sensitization.
  - Monitor the vital signs and check for signs of bleeding.
  - Schedule the patient for an ultrasound in 2 to 4 days to confirm the continued viability of the fetus.
  - Assess the vaginal area for discharge and drainage; note the color and amount.
- PT** Assess and educate the patient for signs of spontaneous abortion (e.g., cramps, bleeding).
- PT** Inform the patient how she can obtain the results from her physician. Make sure that she understands that the results usually will not be available for several weeks. (Results may be available much sooner at major medical centers that perform this test.)
- PT** Educate the patient about how to identify and report signs of endometrial infection (vaginal discharge, fever, crampy abdominal pain).
- PT** Inform the patient about genetic counseling services if needed.

### Abnormal findings

Genetic and biochemical disorders

notes

**chromosome karyotype** (Blood chromosome analysis, Chromosome studies, Cytogenetics, Karyotype)**Type of test** Blood**Normal findings**

Female: 44 autosomes + 2 X chromosomes; karyotype: 46,XX

Male: 44 autosomes + 1 X, 1 Y chromosome; karyotype: 46,XY

**Test explanation and related physiology**

This test is used to study an individual's chromosomal makeup to determine chromosomal defects associated with disease or the risk of developing disease. The term *karyotyping* refers to the arrangement and pairing of cell chromosomes in order from the largest to the smallest to analyze their number and structure. Variations in either can produce numerous developmental abnormalities and diseases. A normal karyotype of chromosomes consists of a pattern of 22 pairs of autosomal chromosomes and a pair of sex chromosomes (XY for the male and XX for the female). Chromosomal karyotype abnormalities can be congenital or acquired. These karyotype abnormalities can occur because of duplication, deletion, translocation, reciprocation, or genetic rearrangement.

Chromosomal karyotyping is useful in evaluating congenital anomalies, mental retardation, growth retardation, delayed puberty, infertility, hypogonadism, primary amenorrhea, ambiguous genitalia, chronic myelogenous leukemia, neoplasm, recurrent miscarriage, prenatal diagnosis of serious congenital diseases (especially in situations of advanced maternal age), Turner syndrome, Klinefelter syndrome, Down syndrome, and other suspected genetic disorders. The products of conception also can be studied to determine the cause of stillbirth or miscarriage.

**Procedure and patient care****Before****PT** Explain the procedure to the patient.

- Determine how the specimen will be collected. Obtain preparation guidelines from the laboratory if indicated.
- Many patients are fearful of the test results and require considerable emotional support.
- In some states, informed consent is required.

**During**

- Specimens for chromosome analysis can be obtained from numerous sources. Leukocytes from a peripheral venipuncture are most easily and most often used for this study.



- Bone marrow biopsies and surgical specimens also can be used sometimes as sources for analysis.
- During pregnancy, specimens can be collected by amniocentesis (see p. 49) and chorionic villus sampling (see p. 241).
- Fetal tissue or products of conception can be studied to determine the reason for loss of a pregnancy.
- Smears and stains of buccal mucosal cells are less costly but not as accurate as other tissue for karyotyping.

### After

- Aftercare depends on how the specimen was collected.
- PT** Inform the patient that test results generally will not be available for several months.
- If an abnormality is identified, often the entire family line must be tested. This can be exhaustive and expensive.
- PT** If the test results show an abnormality, encourage the patient to verbalize his or her feelings. Provide emotional support.

### Abnormal findings

Congenital anomalies  
 Mental retardation  
 Growth retardation  
 Delayed puberty  
 Infertility  
 Hypogonadism  
 Primary amenorrhea  
 Ambiguous genitalia  
 Chronic myelogenous leukemia  
 Neoplasm  
 Recurrent miscarriage  
 Trisomy 21 (Down syndrome)  
 Tay–Sachs disease  
 Sickle cell anemia  
 Turner syndrome  
 Klinefelter syndrome

### notes

***Clostridium difficile*** (*C. diff*, Clostridial toxin assay)**Type of test** Stool**Normal findings**Negative (no *Clostridium* toxin identified)

No other infectious agent identified

**Test explanation and related physiology**

*Clostridium difficile*-associated disease (CDAD) bacterial infections usually affect the intestine (pseudomembranous colitis) and often occur in patients who are immunocompromised or taking broad-spectrum antibiotics (e.g., clindamycin, ampicillin, cephalosporins). The disease can be community acquired, and the severity can range from mild nuisance diarrhea to severe pseudomembranous colitis and bowel perforation.

The infection possibly results from depression of the normal flora of the bowel caused by the administration of antibiotics. The clostridial bacterium produces two toxins (A and B) that cause inflammation and necrosis of the colonic epithelium. *C. difficile* can be diagnosed by obtaining colonic-rectal tissue for this toxin. Stool cultures (see p. 847) for *C. difficile* can be performed but are also labor intensive and take longer to get results.

A rapid detection of *C. difficile* toxin B gene (*tcdB*) in human liquid or soft stool specimens is available. This method rapidly provides a definitive diagnosis of *C. difficile*. Quickly reaching a definitive diagnosis allows CDAD patients to get the proper treatment without delay and reduces hospital stays for inpatients with CDAD. At the same time, they can be placed in isolation sooner to reduce transmission and prevent outbreaks.

A positive PCR result for the presence of the gene-regulating toxin C production (*tcdC*) also indicates the presence of *C. difficile* and toxin A, B, or both.

*Salmonella* spp., *Campylobacter* spp., *Shigella* spp., *Cryptosporidium* spp., *Escherichia coli*, *Vibrio* spp., *Yersinia* spp., and *Giardia lamblia* may also cause acute diarrhea and can be identified in PCR multiplex panel testing. Viral infections caused by noro- or adenoviruses can also be detected on the newer stool testing panels.

**Procedure and patient care****Before**

**PT** Explain the method of stool collection to the patient. Be matter of fact to avoid embarrassment for the patient.

- PT** Instruct the patient not to mix urine or toilet paper with the stool specimen.
- Handle the specimen carefully, as though it were capable of causing infection.

### During

- PT** Instruct the patient to defecate into a clean container. A rectal swab cannot be used because it collects inadequate amounts of stool. The stool cannot be retrieved from the toilet.
- Stool can be obtained from incontinence pads.
  - A stool specimen also can be collected by proctoscopy.
  - Place the specimen in a closed container and transport it to the laboratory to prevent deterioration of the toxin.
  - If the specimen cannot be processed immediately, refrigerate it (depending on laboratory protocol).

### After

- Maintain enteric isolation precautions on all patients until appropriate therapy is completed.

### Abnormal findings

#### ▲ Increased levels

Antibiotic-related pseudomembranous colitis  
*C. difficile* colitis

notes

## coagulating factors concentration (Coagulating factors, Blood-clotting factors)

**Type of test** Blood

### Normal findings

Factor	Normal value (% of “normal”)
II	80-120
V	50-150
VII	65-140
VIII	55-145
IX	60-140
X	45-155
XI	65-135
XII	50-150

### Test explanation and related physiology

These tests measure the quantity of each specific factor thought to be responsible for suspected defects in hemostasis. Testing is available to measure the quantity of the factors listed in the normal findings section. When these factors exist in concentrations below their *minimal hemostatic levels*, clotting is impaired. Deficiencies of these factors may be a result of inherited genetic defects, acquired diseases, or drug therapy. Common medical conditions associated with abnormal factor concentrations are listed in [Table 9](#).

The hemostasis and coagulation system is a homeostatic balance between factors encouraging clotting and those encouraging clot dissolution ([Figure 10](#)). See [Table 10](#) for a list of factor names and routine coagulation test abnormalities associated with factor deficiency.

Fibrinogen (or factor I; see p. 422), like many of the coagulation proteins, is considered an acute reactant protein and is elevated in many severe illnesses. It is also considered a risk factor for coronary heart disease and stroke.

Prothrombin is a vitamin K–dependent clotting factor. Its production in the liver requires vitamin K. This vitamin is fat soluble and is dependent on bile for absorption.

Factor VIII is actually a complex molecule with two components. The first component is related to hemophilia A and is involved in the hemostatic mechanism. The second component is von Willebrand factor and is related to von Willebrand disease. This second component is involved in platelet adhesion and aggregation.

**TABLE 9** Conditions that may result in coagulation factor excess or deficiency

Factor	Increased (excess)	Decreased (deficiency)
I (Fibrinogen)	Acute inflammatory reactions Trauma Coronary heart disease Cigarette smoking	Liver disease (hepatitis or cirrhosis) DIC Congenital deficiency
II (Prothrombin)	ND	Vitamin K deficiency Liver disease Congenital deficiency Warfarin ingestion
V (Proaccelerin)	ND	Liver disease DIC Fibrinolysis
VII (Proconvertin [stable factor])	ND	Congenital deficiency Vitamin K deficiency Liver disease Warfarin ingestion
VIII (Antihemophilic factor)	Acute inflammatory reactions Trauma or stress Pregnancy Birth control pills	Congenital deficiency (e.g., hemophilia A) DIC
von Willebrand factor	ND	Congenital deficiency (e.g., von Willebrand disease) Some myeloproliferative disorders
IX (Christmas factor)	ND	Congenital deficiency (e.g., hemophilia B) Liver disease Nephrotic syndrome Warfarin ingestion DIC Vitamin K deficiency

*Continued*

## 250 coagulating factors concentration

**TABLE 9** Conditions that may result in coagulation factor excess or deficiency—cont'd

Factor	Increased (excess)	Decreased (deficiency)
X (Stuart factor)	ND	Congenital deficiency Liver disease Warfarin ingestion Vitamin K deficiency
XII (Hageman factor)	ND	Congenital deficiency Liver disease DIC

*DIC*, Disseminated intravascular coagulation; *ND*, there is no common disease state known to be associated with an excess of this factor.

*Coagulation factor inhibitors* arise in patients who are congenitally deficient in a specific factor in response to factor replacement therapy. They can also occur spontaneously without known cause or in response to a variety of medical conditions, including the postpartum state, immunologic disorders, certain antibiotic therapies, some malignancies, and old age.

### Interfering factors

- Many of these proteins are heat sensitive, and levels are decreased if the specimen is kept at room temperature.
- Pregnancy or the use of contraceptive medication can increase levels of several of these factors, especially VIII and IX.
- Many of these protein coagulation factors are *acute reactant* proteins. Acute illness, stress, exercise, or inflammation can raise levels.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: blue

### Abnormal findings


See [Table 9](#) on p. 250.

### notes

# Hemostasis and Fibrinolysis

Bold = Vitamin K–dependent coagulating factors


## Reaction 1 (Intrinsic)

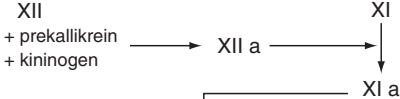


## Reaction 2 (Extrinsic)

**VII + Ca + III (tissue factor)** → Protease


## Reaction 3 (Common pathway)

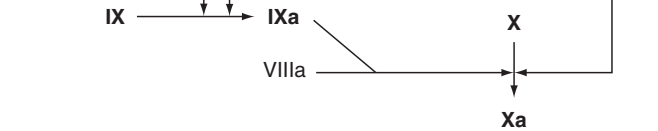


## Reaction 4

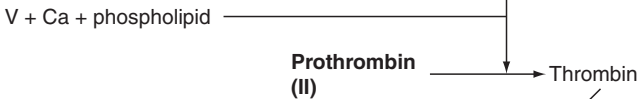

## CLOT STABILIZATION

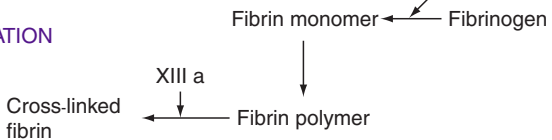

## FIBRINOLYSIS

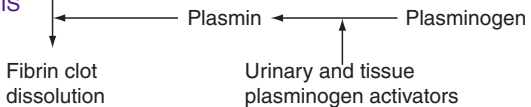

**FIG. 10** Secondary hemostasis (fibrin clot formation) and fibrinolysis (fibrin clot dissolution). Primary hemostasis involves platelet plugging of the injured blood vessel. Secondary hemostasis, as described here, takes place most rapidly on the platelet surface after attachment to the fractured endothelium. Four different reactions result in the formation of fibrin. As seen beneath the dark line in the figure, the fibrin clot supports the platelet clump so that the clot does not get swept away by the tremendous shear forces of the fast-moving blood cells. Fibrinolysis follows formation of the fibrin clot in order to prevent complete occlusion of the injured blood vessel.

**TABLE 10** List of minimum concentration of coagulation factors required for fibrin production<sup>a</sup>

Factor	Name	Quantitation of minimum hemostatic level (mg/dL)	Abnormal coagulation tests associated with deficiency	Blood components to provide specific factor
I	Fibrinogen	60-100	PT, APTT	C, FFP, FWB
II	Prothrombin	10-15	PT	P, WB, FFP, FWB
III	Tissue factor or thromboplastin	Not applicable	PT	
IV	Calcium	See calcium, p. 189		
V	Proaccelerin	5-10	PT, APTT	FFP, FWB
VII	Stable factor	5-20	PT	P, WB, FFP, FWB
VIII	Antihemophilic factor	30	APTT	C, FFP, VIII CONC
IX	Christmas factor	30	APTT	FFP, FWB
X	Stuart factor	8-10	PT, APTT	P, WB, FFP, FWB
XI	Plasma thromboplastin antecedent	25	APTT	P, WB, FFP, FWB
XII	Hageman factor	Yes	APTT	
XIII	Fibrin stabilizing factor	No		P, C, XIII CONC

<sup>a</sup>Recombinant factors are now available for factor VII, VIII, IX, and XIII. Concentrates are also now available for II, VII, VIII, IX, and XIII. *APTT*, Activated partial thromboplastin time; *C*, cryoprecipitate; *FFP*, fresh-frozen plasma; *FWB*, fresh whole blood (< 24 hours old); *P*, unfrozen banked plasma; *PT*, prothrombin time; *VIII CONC*, factor VIII concentrate; *WB*, banked whole blood; *XIII CONC*, factor XIII concentrate.



## cold agglutinins

**Type of test** Blood

### Normal findings

No agglutination in titers  $\leq 1:64$

### Test explanation and related physiology

Cold agglutinins are antibodies (usually IgM) to erythrocytes. All individuals have circulating antibodies directed against RBCs, but their concentrations are often too low to trigger disease (titers  $< 1:64$ ). In individuals with *cold agglutinin disease*, these antibodies are in much higher concentrations. At body temperatures of  $28^{\circ}\text{C}$  to  $31^{\circ}\text{C}$ , such as those encountered during winter months, these antibodies can cause a variety of symptoms, from chronic anemia because of intravascular hemolysis or extravascular sequestration of affected RBCs to acrocyanosis of the ears, fingers, or toes because of local blood stasis in the skin capillaries.

There are two forms of cold agglutinin disease: primary and secondary. The primary form has no precipitating cause. Secondary cold agglutinin disease is a result of an underlying condition, notably *Mycoplasma pneumoniae*. Other possible underlying conditions include influenza, mononucleosis, rheumatoid arthritis, lymphomas, HIV, Epstein–Barr virus, and cytomegalovirus. Temperature regulation is important for the performance of these tests. Under no circumstance should the cold agglutinin specimen be refrigerated.

The cold agglutinin screen is performed on all specimens first to identify most of those with titer values in the normal range. If the result is negative, no titration is required. If the result is positive, a titer with serial saline dilutions is performed.

### Interfering factors

- ✦ Some antibiotics (penicillin and cephalosporins) can interfere with the development of cold agglutinins.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

## 254 cold agglutinins

### Abnormal findings

#### ▲ Increased levels

*Mycoplasma pneumoniae* infection

Viral illness

Infectious mononucleosis

Multiple myeloma

Scleroderma

Cirrhosis

Staphylococemia

Thymic tumor

Influenza

Rheumatoid arthritis

Lymphoma

Systemic lupus erythematosus

Primary cold agglutinin disease

#### notes


## coagulation profile (Coagulation panel)

The coagulation panel is a screening test for abnormal blood clotting. It examines the factors most often associated with a bleeding problem. Abnormal test results can be followed up with other specific tests to confirm or rule out a diagnosis. These tests are listed below and discussed separately. For accuracy, check the normal values at the lab performing the test.

	Normal values
PT (prothrombin time)/INR (international normalized ratio) (see p. 753)	11-12.5 seconds/0.8-1.1
APTT (activated partial thromboplastin time) (see p. 681)	30-40 seconds
Platelet count (see p. 706)	150,000-400,000/mm <sup>3</sup>

This screening profile does not address all causes of bleeding tendencies. Other tests to consider include fibrinogen (see p. 422), D-dimer (see p. 319), coagulation factors (see p. 248), and disseminated intravascular coagulation (DIC) screening (see p. 330).

Test results are often displayed as shown in the illustration below.

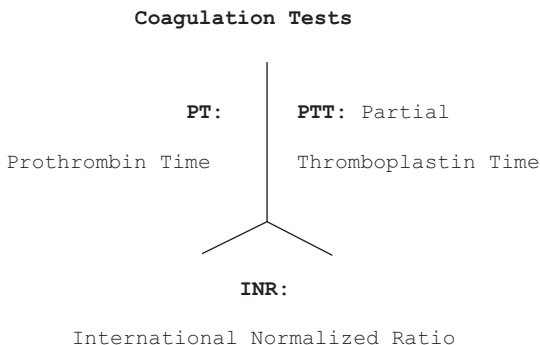

**colon cancer tumor analysis** (Microsatellite instability [MSI] testing, DNA mismatch repair [MMR] genetic testing, BRAF mutation analysis, Oncotype DX colon cancer assay)

**Type of test** Microscopic examination

### Normal findings

Recurrence score < 10

No mismatch repair gene

No microsatellite instability

### Test explanation and related physiology

Patients with stage I colon cancer have a high cure rate with surgery alone. Patients with stage III colon cancer benefit from the use of adjuvant chemotherapy. Patients with stage II colon cancer may or may not benefit from adjuvant chemotherapy. Colon cancer tumor analysis can help differentiate the stage II patients who may benefit from adjuvant chemotherapy. These tests are used to indicate the risk of recurrent colon cancer in the years succeeding surgical treatment.

Deficiencies in DNA mismatch repair (MMR) gene function, because of either decreased gene expression or mutation, result in the accumulation of DNA alterations that can manifest as abnormal shortening or lengthening of microsatellite DNA sequences in the colon cancer cell. This causes microsatellite instability (MSI). Patients with MMR-deficient (MMR-D) colon tumors have high MSI and have been shown to have significantly lower colon cancer recurrence risk. Therefore testing the colon tumor for MMR and MSI can help determine the likelihood of recurrence after surgery and quantify any benefit from adjuvant chemotherapy.

Furthermore, hereditary colon cancers frequently are positive for MSI compared with sporadic colon cancers. Lynch syndrome (a hereditary form of colon cancer) can be suspected if the tumor is MSI positive.

*BRAF* is another important gene that is used to indicate the likelihood that a colon tumor is hereditary. *BRAF* is a kinase-encoding gene in the RAS/RAF/MAPK pathway. The presence of a *BRAF* V600E mutation in a microsatellite unstable tumor indicates that the tumor is probably sporadic and not associated with hereditary nonpolyposis colorectal cancer (HNPCC). The lack of this mutation indicates that a tumor may either be sporadic or HNPCC-associated.

*BRAF* is an important genetic mutation in other cancers such as melanoma, papillary thyroid cancer, hairy-cell leukemia, lung cancer, and other B-cell lymphomas. *BRAF* mutation may be associated with increased risk of recurrence, lymph node metastases, and advanced stage cancer. Vemurafenib is a highly selective and potent inhibitor of *BRAF* V600E. It has marked antitumor effects against melanoma and some other tumors with the *BRAF* V600E mutation. Thus establishing whether *BRAF* mutations exist may be of critical therapeutic importance.

The Oncotype DX colon cancer assay evaluates 12 genes and produces a recurrence score that is closely associated with the patient's risk of recurrent colon cancer 3 years after surgery (the peak time of recurrence). MMR and MSI testing can complement the information provided by the Oncotype DX colon cancer assay.

## Procedure and patient care

### Before

- PT** Inform the patient that an examination for these tumor predictor markers may be performed on their colon cancer tissue.
- PT** Provide psychological and emotional support to the patients colon cancer.

### During

- The surgeon obtains tumor tissue.
- This tissue should be placed on ice or in formalin.
- Part of the tissue is used for routine histology. A portion of the paraffin block is sent to a reference laboratory.

### After

- PT** Explain to the patient that results are usually available in 1 week.

## Abnormal findings

Unfavorable test results indicating a high risk of cancer recurrence

### notes

**colonoscopy, sigmoidoscopy, proctoscopy, anoscopy****Type of test** Endoscopy**Normal findings**

Normal anus, rectum, colon, and distal small bowel

**Test explanation and related physiology**

With fiberoptic colonoscopy, the entire colon from anus to cecum (and often a portion of terminal ileum) can be examined in most patients. *Anoscopy* refers to examination of the anus; *proctoscopy* to examination of the anus and rectum; and *sigmoidoscopy* to examination of the anus, rectum, and sigmoid colon. Sigmoidoscopy can be performed with a rigid ( $\leq 25$  cm from the anus) or flexible ( $\leq 60$  cm from the anus) sigmoidoscope.

Benign and malignant neoplasms, polyps, mucosal inflammation, ulceration, and sites of active hemorrhage can be visualized. Such diseases as cancer, polyps, ulcers, and arteriovenous (AV) malformations also can be visualized. Cancers, polyps, and inflammatory bowel diseases can be biopsied through the scope with cable-activated instruments; sites of active bleeding can be coagulated with the use of laser, electrocoagulation, and injection of sclerosing agents.

This test is recommended for patients who have Hemoccult-positive stools, abnormal sigmoidoscopy, lower GI tract bleeding, abdominal pain, or a change in bowel habits. This test is also recommended for patients who are at high risk for colon cancer. They include patients with a strong personal or family history of colon cancer, polyps, or ulcerative colitis. Colonoscopy is also used for colorectal screening in asymptomatic patients without increased risks for cancer. The U.S. Preventive Services Task Force (USPSTF) recommends screening for colorectal cancer using high-sensitivity fecal occult blood testing, sigmoidoscopy, or colonoscopy beginning at age 50 years and continuing until age 75 years. Virtual colonoscopy is now an option (see p. 273).

Colonoscopies are performed by a physician trained in GI endoscopy in approximately 30 to 60 minutes. They are usually performed in an endoscopy suite or the operating room. Because the patient is sedated, he or she experiences very little discomfort and may not have recall of the procedure. Sigmoidoscopies are performed in about 20 minutes. Sedation is not required. Patients may feel discomfort and the urge to defecate as the scope is inserted.

## Contraindications

- Patients whose medical condition is not stable
- Patients who are bleeding profusely from the rectum. The viewing lens will become covered with blood clots.
- Patients with a suspected perforation of the colon
- Patients with toxic megacolon. These patients may worsen with the test preparation.
- Patients with a recent colon anastomosis (within the past 14-21 days). The anastomosis may break down with significant insufflation of CO<sub>2</sub>.

## Potential complications

- Bowel perforation
- Persistent bleeding from a biopsy site
- Oversedation, resulting in respiratory depression

## Interfering factors

- Poor bowel preparation may result in stool obstructing the lens and precluding adequate visualization of the colon.
- Active bleeding obstructs the lens system and precludes adequate visualization of the colon.

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient.

**PT** Fully inform the patient about the risks of the procedure and obtain an informed consent.

**PT** Assure patients that they will be appropriately draped to avoid unnecessary embarrassment.

**PT** Instruct the patient on appropriate bowel preparation.

A typical preparation for *colonoscopy* may be as follows:

- *7 days before testing*
  - Aspirin or NSAIDS (e.g., Advil, Motrin, Celebrex, Ibuprofen) may be continued. But doctors must provide specific instructions for continuing blood thinners such as Plavix, Pradaxa, clopidogrel, Coumadin, Warfarin, Effient, Prasugrel, or Lovenox.
- *3 days before testing:*
  - Stop eating all nuts, seeds, and popcorn.
- *1 day before testing:*
  - Begin a clear liquid diet. Drink at least 8 glasses of water during the day to avoid dehydration.
  - At noon, take 4 Dulcolax tablets.

## 260 colonoscopy, sigmoidoscopy, proctoscopy, anoscopy

- Mix 64 oz of liquid with 8.3 oz of MiraLAX and place it in the refrigerator.
- At 6 PM, drink one 8-oz glass of the MiraLAX–Gatorade solution and continue drinking one 8-oz glass every 15 minutes thereafter until the mixture is gone.
- PT** Instruct the patient not drink anything colored red, orange, green, or blue because they may interrupt in interpretation and visualization of the intestines.
- *\*A typical preparation for sigmoidoscopy is less severe and may be as follows:*
  - Administer a Fleet enema in the morning of the test.
  - A light breakfast may be ingested.

### During

- Note the following procedural steps for a *colonoscopy*:
  1. IV access is obtained for anesthesia.
  2. After a rectal examination indicates adequate bowel preparation, the patient is sedated.
  3. The patient is placed in the lateral decubitus position, and the colonoscope is placed into the rectum.
  4. Under direct visualization, the colonoscope is directed to the cecum. Often a significant amount of manipulation is required to obtain this position.
  5. As in all endoscopies, air is insufflated to distend the bowel for better visualization.
  6. Complete examination of the large bowel is carried out.
  7. Polypectomy, biopsy, and other endoscopic surgery is performed after appropriate visualization.
  8. When the laser or coagulator is used, the air is removed, and carbon dioxide is used as an insufflating agent to avoid explosion.
- Note the following procedural steps for a *sigmoidoscopy*:
  1. The patient is placed on the endoscopy table or bed in the left lateral decubitus position.
  2. Usually no sedation is required.
  3. The anus is mildly dilated with a well-lubricated finger.
  4. The rigid or flexible sigmoidoscope is placed into the rectum and advanced.
  5. Air is insufflated during the procedure to distend more fully the lower intestinal tract.
  6. The sigmoid, rectum, and anus are visualized.
  7. Biopsies can be obtained and polyps can be removed.



**After**

- PT** Explain to the patient that air has been insufflated into the bowel. She or he may experience flatulence or gas pains.
- Examine the abdomen for evidence of colon perforation (abdominal distention, tenderness, fever, and chills).
  - Monitor vital signs for a decrease in blood pressure and an increase in pulse as an indication of hemorrhage.
  - Inspect the stool for gross blood.
  - Notify the physician if the patient develops increased pain or significant GI bleeding.
  - Allow the patient to eat when fully alert if no evidence of bowel perforation exists.
- PT** Encourage the patient to drink a lot of fluids when intake is allowed. This will make up for the dehydration associated with the bowel preparation.
- PT** Inform the patient that frequent, bloody bowel movements may indicate poor hemostasis after biopsy or polypectomy.
- PT** Educate the patient to report abdominal bloating and inability to pass flatus, which may indicate colon obstruction if a neoplasm was identified.
- PT** Assess the patient for and educate to report weakness and dizziness, which may indicate orthostasis and hypovolemia caused by dehydration.

**Abnormal findings**

Colorectal cancer

Polyps

Inflammatory bowel disease (e.g., ulcerative or Crohn colitis)

Tumors (benign or malignant)

Arteriovenous malformations

Hemorrhoids

Ischemic or postinflammatory strictures

Diverticulosis

**notes**

## colposcopy


**Type of test** Endoscopy

### Normal findings

Normal vagina and cervix

### Test explanation and related physiology

Colposcopy provides an *in situ* macroscopic examination of the vagina and cervix with a *colposcope*, which is a macroscope with a light source and a magnifying lens (Figure 11). With this procedure, tiny areas of dysplasia, carcinoma *in situ* (CIN), and invasive cancer that would be missed by the naked eye can be visualized, and biopsy specimens can be obtained. The study is performed on patients with abnormal vaginal epithelial patterns, cervical lesions, abnormal Pap smears, or positive HPV results

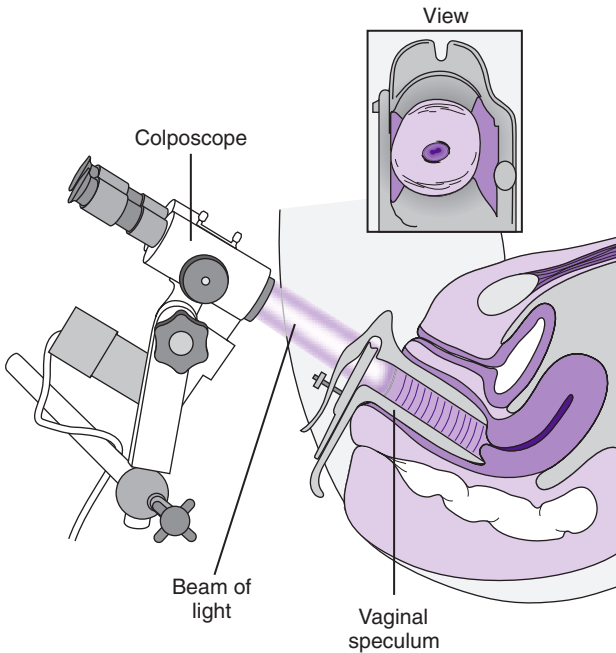

**FIG. 11** Colposcopy. A colposcope is used to evaluate patients with an abnormal Pap smear and a grossly normal cervix.

and on those exposed to diethylstilbestrol *in utero*. This procedure is used to determine the need for cone biopsy (removal and examination of a cone of tissue from the cervix) in evaluating the cause of abnormal cervical cytologic findings.

It is important to realize that colposcopy is useful only in identifying a suspicious lesion. Definitive diagnosis requires biopsy of the tissue. One of the major advantages of this procedure is the capability of directing the biopsy to the area most likely to be truly representative of the lesion. A biopsy performed without colposcopy may not necessarily be representative of the lesion's true pathologic condition, resulting in a significant risk of missing a serious lesion.

The patient will need to have diagnostic conization if

- Colposcopy and endocervical curettage do not explain the problem or match the cytologic findings of the Pap smear within one grade.
- The entire transformation zone is not seen.
- The lesion extends up the cervical canal beyond the vision of the colposcope.

The need for up to 90% of cone biopsies is eliminated by an experienced colposcopist. Endocervical curettage or brushings may accompany colposcopy to detect unknown lesions in the endocervical canal.

There are other promising methods to examine the cervix for neoplasm. *Cervicography* is a procedure in which the cervix is swabbed with an acetic acid solution to identify acetowhite changes in the cervix. A photograph of the cervix is taken with a special camera (Cerviscope) and is sent for imaging analysis for the presence of atypia or metaplasia, intraepithelial neoplasia, or cancer. *Speculoscopy* (PapSure) uses a chemiluminescent light to aid naked-eye or magnified visualization of acetowhite changes on the cervix. *Video colpography* (video colposcopy) has been used for imaging the vagina and cervix and has been proposed for use as a method of cervical cancer screening. In this procedure, a video camera is used to create computerized digital images of the cervix, vaginal fornices, and endocervical canal. The images are then evaluated for signs of cervical cancer. Video colpography may be used in teaching, auditing, and screening of women with low-grade Pap smear abnormalities. *Spectroscopy* optical cervical imaging system is used as an adjunct to colposcopy to identify areas of the cervix with the highest likelihood of high-grade CIN on biopsy. Spectroscopy shines a light on the cervix and analyzes how different areas of the cervix respond to the light. The system produces a color map that distinguishes between healthy

## 264 colposcopy

and potentially diseased tissue to indicate where biopsy samples should be taken.

Finally, there are newer methods of gene identification using methylation markers and qualitative fluorescence *in situ* hybridization (FISH) methods to identify abnormal genes associated with early markers of neoplastic cervical disease.

### Contraindications

- Patients with heavy menstrual flow

### Interfering factors

- Failure to cleanse the cervix of foreign materials (e.g., creams, medications) may impair visualization.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Obtain informed consent if required by the institution.

#### During

- Note the following procedural steps:
  1. The patient is placed in the lithotomy position, and a vaginal speculum is used to expose the vagina and cervix.
  2. The cervix is cleansed with a 3% acetic acid solution to remove excess mucus and cellular debris. The acetic acid also accentuates the difference between normal and abnormal epithelial tissues. Abnormal epithelium becomes white with application of dilute acetic acid.
  3. An aggressive ecto-endocervical Pap smear using curettage is then performed.
  4. The colposcope is focused on the cervix, which is then carefully examined.
  5. Usually the entire lesion can be outlined and the most atypical areas selected for biopsy specimen removal.
- Note that colposcopy is performed by a physician, nurse practitioner, or physician's assistant in approximately 5 to 10 minutes.
- PT** Tell the patient that some women complain of pressure pains from the vaginal speculum and that momentary discomfort may be felt if biopsy specimens are obtained.

#### After

**PT** Inform the patient that she may have vaginal bleeding if biopsy specimens were taken. Suggest that she wear a sanitary pad.

- PT** Instruct the patient to abstain from intercourse and to not insert anything (except a tampon) into the vagina until healing of a biopsy is confirmed.
- PT** Inform the patient about when and how to obtain results of this study.

### **Abnormal findings**

Dysplasia

Carcinoma *in situ*

Invasive cancer

Cervical lesions

**notes**

**complement assay** (Total, C3 and C4 complement)**Type of test** Blood**Normal findings**

Total complement: 30-75 units/mL

C3: 75-175 mg/dL

C4: 22-45 units/mL

**Test explanation and related physiology**

Measurements of complement are used primarily to diagnose hereditary and acquired deficiencies of complement peptides and to monitor the activity of infectious or autoimmune diseases (e.g., systemic lupus erythematosus, nephritis, membranoproliferative nephritis, poststreptococcal nephritis).

Serum complement is a group of 31 proteins that act as enzymes, cofactors, inhibitors, and membrane-integrated proteins. These effect a cascade-like series of reactions that lead to the synthesis of a group of proteins that facilitate the immunologic and inflammatory responses. The total complement, sometimes labeled *CH 50*, is made up of a series of reactions involving proteins C1 through C9 (classic cascade reactions). Besides these major components, subcomponents are involved in the system. When activated, total complement (and some precursor proteins) acts to increase vascular permeability, allowing antibodies and WBCs to be delivered to the area of the immune/antigen complex. Complement also acts to increase chemotaxis (attracting WBCs to the area), phagocytosis, and immune adherence of the antibody to antigen. These processes are vitally important in the normal inflammatory or immunologic response.

Reduced complement levels can be congenital or acquired. As the complement system is activated, the complement components are consumed or used up. If the system is persistently or overly activated, serum levels can fall. The complement system is instigated by the presence of antibody-antigen complexes. As in hereditary angioedema, complement components are used up, and serum levels fall. Diseases associated with these immune complexes include serum sickness, lupus erythematosus, infectious endocarditis, renal transplant rejection, vasculitis, some forms of glomerulonephritis, and infections. As these diseases are successfully treated, complement levels can be expected to return to normal.

The total complement assay should be used as a screen for suspected complement-related diseases before ordering individual

complement component assays. A deficiency of an individual component of the complement cascade will result in an undetectable total complement level. For a list of common diseases associated with complement abnormalities, see [Table 11](#). Note, however, that this list is not complete. Complement abnormalities may occur in the face of normal blood levels when particular complement proteins are not functioning properly.

Complement levels can also be measured in other bodily fluids such as pleural, pericardial, and synovial fluids. Low fluid complement levels are characteristic of effusions from patients with rheumatoid arthritis (despite elevated serum levels), systemic lupus erythematosus, and bacterial infections.

### Interfering factors

- C3 is very unstable at room temperature.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**TABLE 11** Diseases associated with complement deficiencies

Complement deficiency	Associated disease
C1q	Recurrent bacterial infection
C1r	Discoid lupus, glomerulonephritis
C1s	Systemic lupus
C1-INH	Hereditary angioedema
C1	Autoimmune diseases, hypogammaglobulinemia
C2	Lupus, glomerulonephritis, recurrent bacterial infections
C3	Recurrent bacterial infections
C4	Systemic lupus
C5	Systemic lupus, recurrent infections, <i>Neisseria</i> infection
C6	<i>Neisseria</i> infections
C7	Scleroderma, <i>Neisseria</i> infections, rheumatoid arthritis
C8	<i>Neisseria</i> infections
C9	<i>Neisseria</i> infections

## Abnormal findings

### ▲ Increased levels

- Rheumatic fever (acute)
- Myocardial infarction (acute)
- Ulcerative colitis
- Cancer

### ▼ Decreased levels

- Cirrhosis
- Autoimmune disease (e.g., systemic lupus erythematosus)
- Serum sickness (immune complex disease)
- Glomerulonephritis
- Lupus nephritis
- Renal transplant rejection (acute)
- Protein malnutrition
- Anemia
- Malnutrition
- Hepatitis
- Rheumatoid arthritis
- Sjögren syndrome
- Severe sepsis

notes



## complete blood count and differential count (CBC and diff)

The CBC and differential count are a series of tests of the peripheral blood; they provide a tremendous amount of information about the hematologic system and many other organ systems. These tests are inexpensively, easily, and rapidly performed as a screening test. The CBC and differential count include automated measurement of the following studies, which are discussed separately:

RBC count (see p. 770)

Hemoglobin (Hgb; see p. 488)

Hematocrit (Hct; see p. 485)

Blood smear (see p. 149)

Platelet count (see p. 706)

Platelet volume, mean (MPV) (see p. 712)

RBC indices (see p. 772)

Mean corpuscular volume (MCV)

Mean corpuscular hemoglobin (MCH)

Mean corpuscular hemoglobin concentration (MCHC)

Red blood cell distribution width (RDW)

WBC count and differential count (see p. 974)

Neutrophils (polymorphonuclear leukocytes or *polys*, segmented cells or *segs*, band cells, stab cells)

Lymphocytes

Monocytes

Eosinophils

Basophils

### notes

**comprehensive metabolic panel** (CMP, Chem 12, Chemistry panel, SMA 12 (Sequential Multiple Analysis), SMA 20)

The CMP is a broad screening tool used to evaluate organ function and to check for conditions such as liver disease, kidney disease, and diabetes. It is valuable for monitoring known conditions, such as hypertension, and for monitoring kidney- and liver-related side effects of medications. The CMP is often part of the blood workup for medical examinations or yearly physical examinations. Depending on the reason for the test, blood may be drawn after fasting or on a random basis. Abnormal test results can be followed up with other specific tests to rule out or confirm a suspected diagnosis. These tests are listed below and are discussed separately. Sample normal values are shown. For accuracy, use ranges from the lab performing the test.

	Normal values
<b>Electrolytes</b>	
Sodium (see p. 835)	136-145 mEq/L
Potassium (see p. 724)	3.5-5 mEq/L
Chloride (see p. 233)	98-106 mEq/L
CO <sub>2</sub> (carbon dioxide, HCO <sub>3</sub> <sup>-</sup> ; see p. 197)	23-30 mEq/L
Calcium, total (see p. 189)	9-10.5 mg/dL
<b>Glucose</b>	
Glucose (see p. 462)	74-106 mg/dL
<b>Kidney tests</b>	
BUN (blood urea nitrogen; see p. 155)	10-20 mg/dL
Creatine (see p. 301)	0.5-1.1 mg/dL
<b>Liver tests</b>	
ALP (alkaline phosphatase; see p. 29)	30-120 U/L
ALT (alanine aminotransferase; see p. 21)	4-36 U/L
AST (aspartate aminotransferase; see p. 125)	0-35 U/L
Albumin (see p. 746)	3.5-5 g/dL
Total protein (see p. 746)	6.4-8.3 g/dL
Bilirubin (see p. 137)	
Total	0.3-1 mg/dL
Direct	0.1-0.3 mg/dL

Test results are often displayed as illustrated below.

Na	Cl	BUN	Gluc
K	CO <sub>2</sub>	Creat	
Ca	TP	AST	Bili
PO <sub>4</sub>	Alb	ALT	

**computed tomography of the abdomen and pelvis** (CT scan of the abdomen and pelvis, Helical/spiral CT scan of the abdomen and pelvis)

**Type of test** X-ray with or without contrast dye

**Normal findings**

No evidence of abnormality

**Test explanation and related physiology**

The CT scan of the abdomen is a noninvasive, yet very accurate, x-ray procedure used to diagnose pathologic conditions, such as tumors, cysts, abscesses, inflammation, perforation, bleeding, obstruction, aneurysms, and calculi of the abdominal and retroperitoneal organs. Repeating the CT scan after IV and/or oral administration of an iodine-containing contrast dye can enhance the CT image.

Liver tumors, abscesses, trauma, cysts, and anatomic abnormalities can be seen; pancreatic tumors, pseudocysts, inflammation, calcification, bleeding, and trauma can also be detected. The kidneys and urinary outflow tract are well visualized.

Renal tumors and cysts, ureteral obstruction, stones, and congenital renal and ureteral abnormalities are easily seen with the use of IV contrast injection. Calculi can be seen without IV contrast. Extravasation of urine secondary to trauma or obstruction can also be demonstrated easily. Adrenal tumors and hyperplasia are best diagnosed with this technique.

Large tumors, perforations of the bowel, and appendicitis can be identified with the CT scan, especially when oral contrast is ingested. The spleen can be well visualized for hematoma, laceration, fracture, tumor infiltration, and splenic vein thrombosis with CT scanning. The retroperitoneal lymph nodes can be evaluated. These are usually present, but all nodes with a diameter greater than 2 cm are considered abnormal. The abdominal aorta and its major branches can be evaluated for aneurysmal dilation and intramural thrombi. The pelvic structures (including the uterus, ovaries, fallopian tubes, prostate, and rectum) and musculature can be evaluated for tumors, abscesses, infection, or hypertrophy. Ascites and hemoperitoneum can easily be demonstrated with the CT scan. Tumors, abscesses, or perforation of the pelvic organs can be seen when the CT scan is directed to the pelvis. Perineal CT scanning can demonstrate perianal abscesses or perirectal tumors or infection.

*Helical* (also called *spiral*, or *volume-averaging*) CT scanning can obtain many images in a very short period of time. The entire abdomen can be scanned in slightly more than a few seconds with one breath hold. The slices are very thin (1-5 mm). With thin slices and rapid accession, breathing and motion distortion are minimized. This produces faster and more accurate images.

With this CT technique, 200 to 500 individual images can be obtained. Volume imaging with three-dimensional (3D) real-time display of the data allows the interpreter to visualize and analyze the data in three dimensions. With these advances in software, two-dimensional (2D) and 3D reconstructions of data can provide very accurate images of the intraabdominal organs and especially the mesenteric vessels in a few seconds. This allows radiologists to see these structures from multiple directions.

A 3D perspective can be added to the abdominal and pelvic organs or tumors that are imaged. This provides data for virtual colonoscopy and virtual angiography.

*Virtual colonoscopy* (or *CT colonography*) uses a CT scanner and virtual reality software to look inside the body without inserting a colonoscope (for conventional colonoscopy see p. 258). Virtual colonoscopy is an appropriate alternative to a screening endoscopic colonoscopy. No sedation is required, and no discomfort is experienced. Patients need a cleansing bowel preparation before the test. The virtual colonoscopy procedure takes place in the radiology department. It begins with the insertion of a small, flexible rubber tube in the rectum. Air is inserted through this tube to inflate the colon for better visualization. The air acts as a contrast medium. The test is completed in 10 to 20 minutes. Because no sedation is required, patients are free to leave the CT suite without the need for observation and recovery. Patients can resume normal activities after the procedure, and they can eat, work, or drive without a delay. Unlike endoscopic colonoscopy, polypectomy or biopsy cannot be performed with virtual testing. If abnormalities are found with virtual colonoscopy, conventional colonoscopy is needed.

Increasingly, *fusion CT/positron emission tomography (PET)* scans are now being used to provide both anatomic and physiologic information that can be fused onto one image. This allows the image to locate pathology and to indicate whether it is benign or malignant. Fusion CT/PET scans can provide an accurate image of the entire colon and indicate whether any abnormality that is seen is malignant.

*CT arteriography or angiography* is performed through the use of multichannel helical CT scanning. After IV injection of

contrast, CT imaging can demonstrate the arteries in any given organ. With computerized subtraction of the surrounding tissue, the arteries can be even better displayed. 3D recreations of the aorta and other abdominal vessels are possible. This is particularly helpful in identifying renal artery stenosis and hepatic vasculature for cancer-related resections. Renal helical arteriography can be used to demonstrate and evaluate each functional phase of urinary excretion. This procedure is also called *dynamic CT scanning*.

*CT nephrotomography* can be done by computerized recreation of a 3D image of the kidneys, renal pelvis, and ureters. This is particularly helpful in identifying ureteral stones and small tumors of the kidney or collecting system. Using different protocols and radiopaque contrast, kidney function can be evaluated. This does require significant radiation exposure. A different protocol designed to identify ureteral stones can be performed with very little radiation exposure. This is called *CT urography*.

With the increasing use and development of 3D volumetric imaging, radiologists have expanded CT scanning to assist pathologists, coroners, and medical examiners to investigate cadavers for clues as to the cause of death. This is now being termed *virtual autopsy*. This includes CT or MRI whole-body, postmortem imaging. With these techniques, image-directed biopsies can be performed to obtain tissue for the pathologists to review. Postmortem angiograms can be performed to more accurately indicate occlusive disease that may have contributed to death.

The CT scan can be used as a guide to aspirate fluid from the abdomen or one of the abdominal organs. This fluid can be sent for cultures and other studies. The CT scan can also be used to guide biopsy needles into areas of abdominal tumors to obtain tissue for study. Catheters for drainage of intraabdominal abscess can be placed with CT guidance.

The CT scan is an important part of staging and monitoring of many tumors before and after therapy. Tumors of the colon, rectum, liver, breast, lung, prostate, ovary, uterus, kidney, lymph, and adrenal gland commonly recur in the abdomen. Recurrence can be detected early with the CT scan.

### **Contraindications**

- Patients who are allergic to iodinated dye
- Patients who are pregnant, unless the benefits outweigh the risks
- Patients whose vital signs are unstable
- Patients who are very obese (usually > 400 lb)
- Patients who are claustrophobic

## Potential complications

- Allergic reaction to iodinated dye
- See p. xxi for appropriate interventions concerning the care of patients with iodine allergy.
- Acute renal failure from dye infusion. Adequate hydration before the infusion may reduce this likelihood.
- Lactic acidosis may occur in patients who are taking metformin and receiving iodine contrast. Metformin should be held the day of the test to prevent this complication.

## Interfering factors

- Presence of metallic objects (e.g., hemostasis clips)
- Retained barium from previous studies
- Large amounts of fecal material or gas in the bowel

## Procedure and patient care

### Before

- PT** Explain the procedure. The patient's cooperation is necessary because he or she must lie still during the procedure. See p. xviii for radiation exposure and risks.
- Obtain informed consent if required by the institution.
  - Assess the patient for allergies to iodinated dye.
- PT** Show the patient a picture of the CT machine. Encourage the patient to verbalize his or her concerns because some patients may have claustrophobia. Most patients who are mildly claustrophobic can be scanned after appropriate premedication with anti-anxiety drugs.
- Keep the patient NPO (nothing by mouth) for at least 4 hours before the test if oral contrast is to be administered; however, this test can be performed on an emergency basis on patients who have recently eaten.
  - Note that this procedure is usually performed by a radiologist in less than 30 minutes. If dye is administered, the procedure time may be doubled because the CT scan is done both with and without contrast dye.
- PT** Tell the patient that the discomforts associated with this study include lying still on a hard table and the peripheral venipuncture. Mild nausea is a common sensation when contrast dye is used. An emesis basin should be readily available. Some patients may experience a salty taste, flushing, and warmth during the dye injection.

## 276 computed tomography of the abdomen and pelvis

### During

- Note the following procedural steps for an abdominal CT scan:
  1. The patient is taken to the radiology department and placed on the CT scan table.
  2. The patient then is placed in an encircling body scanner (gantry). The x-ray tube travels around the gantry and takes pictures of the various levels of the abdomen and pelvis. Any motion causes blurring and streaking of the final picture. Therefore the patient is asked to remain motionless during x-ray exposure. Computer monitoring equipment allows for immediate display of the CT scan image, which is then recorded digitally. In a separate room, the technicians manipulate the CT scan table to change the level of the abdomen that is scanned.
  3. Through audio communication, the patient is instructed to hold his or her breath during x-ray exposure.
- Remember that oral and IV iodinated x-ray contrast dye provides better results for this test. One can accurately differentiate the GI organs from the other abdominal organs with oral contrast. Likewise, the vessels and ureters are contrasted with the surrounding structures with use of IV dye. Sometimes contrast is given rectally to visualize the pelvic organs.

### After

- PT** Encourage the patient to drink fluids to avoid dye-induced renal failure and to promote dye excretion.
- PT** Inform the patient that diarrhea may occur after ingestion of the oral contrast.
- Evaluate the patient for delayed reaction to dye.



## Abnormal findings

### Liver

Tumor  
Abscess  
Bile duct dilation

### Pancreas

Tumor  
Pseudocyst  
Inflammation  
Bleeding

### Spleen

Hematoma  
Fracture  
Laceration  
Tumor  
Venous thrombosis

### Gallbladder/biliary system

Gallstones  
Tumor  
Bile duct dilation

### Kidneys

Tumor  
Cyst  
Ureteral obstruction  
Calculi  
Congenital abnormalities  
Renal artery stenosis

### Adrenal gland

Adenoma  
Cancer  
Pheochromocytoma  
Hemorrhage  
Myelolipoma hyperplasia

### GI tract

Perforation  
Tumor  
Inflammatory bowel disease  
Diverticulitis  
Appendicitis

### Uterus, tubes, ovaries

Tumor  
Abscess  
Infection  
Hydrosalpinx  
Cyst  
Fibroid

### Prostate

Hypertrophy  
Tumor

### Retroperitoneum

Tumor  
Lymphadenopathy

### Other

Abdominal aneurysm  
Ascites, hemoperitoneum  
Abscess

notes

## computed tomography of the brain (CT scan of the brain, Helical/spiral CT scan of the brain)

**Type of test** X-ray with or without contrast dye

### Normal findings

No evidence of pathologic conditions

### Test explanation and related physiology

CT of the brain consists of a computerized analysis of multiple tomographic x-rays taken of the brain tissue at successive layers, providing a 3D view of the cranial contents. The CT image provides a view of the head as if one were looking down through its top. The variation in density of each tissue allows for variable penetration of the x-ray beam. An attached computer calculates the amount of x-ray penetration of each tissue and displays this as shades of gray. This is then displayed digitally on a computer monitor as a series of actual anatomic pictures of coronal sections of the brain.

The CT scan is used in the differential diagnosis of intracranial neoplasms, cerebral infarctions, ventricular displacement or enlargement, cortical atrophy, cerebral aneurysms, intracranial hemorrhage and hematoma, and AV malformation. Information about the ventricular system can be obtained by CT scanning. Multiple sclerosis and other degenerative abnormalities can be identified also.

Visualization of a neoplasm, previous infarction, or any pathologic process that destroys the blood-brain barrier may be enhanced by IV injection of an iodinated contrast dye. CT scans may be repeated frequently to monitor the progress of any disease or to monitor the healing process. In most cases, CT scanning has eliminated the need for more invasive procedures, such as *cerebral arteriography* and *pneumoencephalography*. MRI brain scanning can provide more and different information in some instances and may be used in place of CT scanning of the brain.

*CT arteriography* is performed immediately after venous contrast injection. 3D recreations of the carotid artery and its branches are also extremely helpful in the evaluation of cerebral vascular disease.

*Helical* (also called *spiral*, or *volume-averaging*) *CT scanning* produces faster and more accurate images. Because the helical CT scan images the selected area in less than

30 seconds, the entire study can be performed with one breath hold. Therefore breathing and motion misrepresentations are eliminated. This is particularly helpful in scanning uncooperative adults or children. Through volume averaging, 3D images can be recreated. Furthermore, when contrast material is used, contrast imaging is markedly improved.

### Contraindications

- Patients who are allergic to iodinated dye
- Patients who are claustrophobic
- Patients who are pregnant, unless benefits outweigh risks
- Patients whose vital signs are unstable
- Patients who are very obese (usually > 400 lb)

### Potential complications

- Allergic reaction to iodinated dye: See p. xxi for appropriate interventions concerning the care of patients with iodine allergy.
- Acute renal failure from dye infusion: Adequate hydration beforehand may reduce this likelihood.
- Lactic acidosis may occur in patients who are taking metformin and receiving iodine contrast. Metformin should be held the day of the test to prevent this complication.
- Apnea if xenon is used because it is an anesthetic gas

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. The patient's cooperation is necessary, because he or she must lie still during the procedure. See p. xviii for radiation exposure and risks.
- Obtain informed consent if required by the institution.
- PT** Show the patient a picture of the CT machine and encourage the patient to verbalize his or her concerns because some patients may have claustrophobia. Most patients who are mildly claustrophobic can be scanned after appropriate premedication with antianxiety drugs.
- Keep the patient NPO for 4 hours before the study if oral contrast is to be used.
- PT** Instruct the patient that wigs, hairpins, clips, or partial dental plates cannot be worn during the procedure because they hamper visualization of the brain.
- Assess the patient for allergies to iodinated dye.
- PT** Tell the patient that he or she may hear a clicking noise as the scanning machine moves around the head.

## 280 computed tomography of the brain

### During

- Note the following procedural steps for the brain CT scan:
  1. The patient lies in the supine position on an examining table with the head resting on a platform. The face is not covered, and the patient can see out of the machine at all times. Sponges are placed along the side of the head to ensure that the patient's head does not move during the study.
  2. The scanner passes an x-ray beam through the brain from multiple angles.
  3. If an iodinated dye will be used, a peripheral IV line is started and the iodinated dye is administered through it. The entire scanning process is repeated.
- Note that this procedure is performed by a radiologist in less than 1 hour. If dye is administered, the procedure time is doubled because the CT scan is done with and without the contrast dye.
- PT** Tell the patient that discomfort associated with this study includes lying still on a hard table and peripheral venipuncture. Mild nausea is a common sensation when contrast dye is used. An emesis basin should be readily available. Some patients may experience a salty taste, flushing, and warmth during the dye injection.

### After

- PT** Encourage the patient to drink fluids because dye is excreted by the kidneys and causes diuresis.

### Abnormal findings

Intracranial neoplasm  
Cerebral infarction  
Ventricular displacement  
Ventricular enlargement  
Cortical atrophy  
Cerebral aneurysm  
Intracranial hemorrhage  
Hematoma  
Arteriovenous malformation  
Meningioma  
Multiple sclerosis  
Hydrocephalus  
Abscess

### notes

## computed tomography of the chest (Chest CT scan, Helical/spiral CT scan of the chest)

**Type of test** X-ray with or without contrast dye

### Normal findings

No evidence of pathologic conditions

### Test explanation and related physiology

CT of the chest is a noninvasive yet very accurate x-ray procedure for diagnosing and evaluating pathologic conditions such as tumors, nodules, hematomas, parenchymal coin lesions, cysts, abscesses, pleural effusion, and enlarged lymph nodes affecting the lungs and mediastinum. Tumors, cysts, and fractures of the chest wall and pleura can also be seen. When IV contrast material is given, vascular structures can be identified, and a diagnosis of aortic or other vascular abnormality can be made. With oral contrast, the esophagus and upper structures can be evaluated for tumor and other conditions. This procedure provides a cross-sectional view of the chest and is especially useful in detecting small differences in tissue densities, demonstrating lesions that cannot be seen with conventional radiology and tomography.

The entire chest can be scanned in slightly more than a few seconds with one breath hold. The slices are very thin (1-5 mm). With thin slices and rapid acquisition, breathing and motion distortion are minimized. This produces faster and more accurate images.

With this CT study, 200 to 500 individual images can be obtained. Volume imaging with 3D real-time display of the data allows the interpreter to visualize and analyze the data in three dimensions. With these advances in software, 2D and 3D reconstructions of data can provide very accurate images of the heart (see p. 284), lungs, chest wall, pleura, esophagus, great vessels, and soft tissue in a few seconds, allowing radiologists to see these structures from multiple directions. With this new technology, *virtual bronchoscopy* and *virtual esophagoscopy* may increasingly be used in place of their invasive counterparts.

*Four-dimensional computed tomography (4DCT)* has also been developed using 3D imaging and adding time (fourth dimension) sequences in relation to respiratory motion or time from injection of contrast media. This is particularly helpful in studying a moving structure, such as the heart. Breathing no longer introduces artifactual errors in imaging of the organs of the chest (or abdomen). 4DCT is especially helpful in *CT scanning of the*

## 282 computed tomography of the chest

*parathyroid glands.* In this scan, 3D images are obtained at different times (fourth dimension) after injection of contrast. This has been extremely beneficial for patients whose parathyroid glands are not visualized by ultrasonography or parathyroid scan (see p. 677). This technique is also used to plan preoperative parathyroid surgery.

### Contraindications

- Patients who are pregnant unless the benefits outweigh the risks
- Patients who are allergic to iodinated dye
- Patients who are claustrophobic
- Patients who are very obese (usually > 400 lb)
- Patients whose vital signs are unstable

### Potential complications

- Allergic reaction to iodinated dye. See p. xxi for appropriate interventions concerning the care of patients with iodine allergy.
- Acute renal failure from dye infusion. Adequate hydration beforehand may reduce this likelihood.
- Lactic acidosis may occur in patients who are taking metformin and receiving iodine contrast. Metformin should be held the day of the test to prevent this complication.

### Procedure and patient care

#### Before

**PT** Explain the procedure. The patient's cooperation is necessary because he or she must lie still during the procedure. See p. xviii for radiation exposure and risks.

- Obtain informed consent if required by the institution.
- Assess the patient for allergies to iodinated dye.

**PT** Show the patient a picture of the CT machine and encourage the patient to verbalize concerns regarding claustrophobia. Most patients who are mildly claustrophobic can tolerate this study after appropriate premedication with antianxiety drugs.

- Keep the patient NPO for 4 hours before the test in the event that contrast dye is administered.

#### During

- Note the following procedural steps for the chest CT scan:
  1. The patient is taken to the radiology department and asked to remain motionless in a supine position because any motion will cause blurring and streaking of the final picture.

2. An encircling x-ray camera (body scanner) takes pictures at varying intervals and levels over the chest area. Often IV dye is administered to enhance the chest image, and the x-ray studies are repeated.
- Note that this procedure is performed by a radiologist in less than 30 minutes.
- PT** Tell the patient that the discomforts associated with this study include lying still on a hard table and peripheral venipuncture. Mild nausea is a common sensation when contrast dye is used. An emesis basin should be readily available. Some patients may experience a salty taste, flushing, and warmth during the dye injection.

### After

- PT** Encourage patients who received dye injection to increase their fluid intake because the dye is excreted by the kidneys and causes diuresis.

### Abnormal findings

Pulmonary tumor

Inflammatory nodules

Granuloma

Cyst

Pulmonary emboli

Pleural effusion

Enlarged lymph nodes

Aortic aneurysm

Postpneumonic scanning

Pneumonitis

Esophageal tumor

Hiatal hernia

Mediastinal tumor (e.g., lymphoma, thymoma)

Primary or metastatic chest wall tumor

### notes

**computed tomography of the heart** (Coronary CT angiography, Coronary calcium score)

**Type of test** X-ray with contrast dye

**Normal findings**

No evidence of coronary stenosis; calcium score average for age and gender

**Test explanation and related physiology**

This test is being used to help stratify patients according to risks of future cardiac events, to instigate preventive medicinal interventions (e.g., statin drugs), to monitor progression of coronary vascular disease and effects of statin drugs, to evaluate chest pain, and to indicate the need for stress testing or coronary angiography.

The entire heart can be scanned in 10 seconds with one breath hold. The slices are very thin (1-5 mm). With thin slices and rapid accession, breathing and motion distortion are minimized. With advances in software technology, 2D and 3D reconstructions of data can provide very accurate images of the heart and coronary vessels in a few seconds, allowing radiologists to see these structures from multiple directions. Furthermore, with shorter scanning times, IV contrast effect can be greater while using less contrast volume. The scanner can allow routine cardiac gating that synchronizes the scanning with each heartbeat, thereby eliminating further motion distortion.

Cardiac CT provides evaluation of the myocardium, cardiac chambers, cardiac valves, coronary arteries, and great vessels and for detection of pulmonary emboli. Calcified atheromatous plaques can be seen and quantified (calcium score) with this study. Coronary calcium is a surrogate marker for coronary atherosclerotic plaque. In the coronary arteries, calcifications occur almost exclusively in the context of atherosclerotic changes. Within a coronary vessel or larger segment of the vessel, the amount of coronary calcium correlates moderately closely with the extent of atherosclerotic plaque burden.

The *Agatston score* has most frequently been used to quantify the amount of coronary calcium in CT. In the majority of asymptomatic men older than 55 years of age and women older than 65 years of age, calcification can be detected. See [Table 12](#) for categorizing absolute Agatston scores. It is well established that individuals with Agatston scores above 400 have an increased occurrence of coronary procedures (bypass, stent placement,



**TABLE 12** Agatston score categories for coronary calcium quantification

	Minimal	Moderate	Increased	Severe
Agatston score	< 10	11-99	100-400	> 400

angioplasty) and events (myocardial infarction and cardiac death) within 2 to 5 years after the test.

This test can directly and accurately visualize coronary artery lumen after IV injection of a contrast agent (*coronary CT angiography*). Regular and low heart rates are a prerequisite for reliable visualization of the coronary arteries.

### Contraindications

- Patients who are pregnant
- Patients who are allergic to iodinated dye
- Patients who are very obese (usually > 400 lb)
- Patients whose vital signs are unstable

### Potential complications

- Allergic reaction to iodinated dye. See p. xxi for appropriate interventions concerning the care of patients with iodine allergies.
- Acute renal failure from dye infusion. Adequate hydration beforehand may reduce this likelihood.
- Lactic acidosis may occur in patients who are taking metformin and receiving iodine contrast. Metformin should be held the day of the test to prevent this complication.

### Procedure and patient care

#### Before

- PT** Explain the procedure. The patient's cooperation is necessary because he or she must lie still during the procedure. See p. xviii for radiation exposure and risks.
- Obtain informed consent if required by the institution.
  - Assess the patient for allergies to iodinated dye.
  - Assess the patient's vital signs. If the heart rate exceeds protocol levels, administer a rapid-acting beta-blocker or angiotensin-converting enzyme (ACE) inhibitor per protocol orders.
- PT** Show the patient a picture of the CT machine and encourage the patient to verbalize concerns regarding claustrophobia.
- Keep the patient NPO for 4 hours before the test.

## 286 computed tomography of the heart

### During

- Note the following procedural steps for the cardiac CT scan:
  1. The patient is taken to the CT department and asked to remain motionless in a supine position.
  2. ECG leads are applied to synchronize the ECG signal to the image data (gating).
  3. An encircling x-ray camera (body scanner) takes pictures at varying intervals and levels over the heart while the patient holds his or her breath (for about 10 seconds).
  4. A nonenhanced scan is performed first for calcium scoring.
  5. If the calcium scoring is below threshold levels of the protocol, IV dye is rapidly administered through a large-bore IV catheter and the scan is repeated.
  6. A fast-acting nitrate (usually nitroglycerin) is administered to maximize coronary dilation.
- Note that a radiologist or cardiologist performs this procedure in about 20 minutes.
- PT** Tell the patient that discomfort associated with this study includes lying still on a hard table and peripheral venipuncture. Nausea is a common sensation when contrast dye is used. An emesis basin should be readily available. Some patients may experience a salty taste, flushing, and warmth during the dye injection.

### After

- PT** Encourage patients to increase their fluid intake because the dye is excreted by the kidneys and causes diuresis.
- PT** Tell the patient that a headache from the nitroglycerin is not uncommon.

### Abnormal findings

Coronary vascular disease  
Coronary vascular congenital anomalies  
Ventricular aneurysm  
Aortic aneurysm or dissection  
Pulmonary emboli  
Cardiac tumors  
Myocardial scarring  
Cardiac valvular disease

### notes

## Coombs test, direct (Direct antiglobulin test [DAT])

**Type of test** Blood

### Normal findings

Negative; no agglutination

### Test explanation and related physiology

This test is performed to identify immune hemolysis (lysis of RBCs) or to investigate hemolytic transfusion reactions. Most of the antibodies to RBCs are directed against the ABO/Rh blood grouping antigens, such as those that occur in hemolytic anemia of the newborn or transfusion of incompatible blood. When a transfusion reaction occurs, the Coombs test can detect the patient's antibodies or complement components coating the transfused RBCs. Therefore the Coombs test is very helpful in evaluating suspected transfusion reactions.

Non-blood grouping antigens can develop on the RBC membrane and stimulate formation of antibodies. Such drugs as levodopa or penicillin cause this. Also, in autoimmune diseases, antibodies not originally directed against the patient's RBCs can attach to the RBCs and cause hemolysis that is detected by the direct Coombs test. Frequently, the inciting factor for the production of these autoantibodies against RBCs is not associated with any identifiable disease, and the resulting hemolytic anemia is called *idiopathic*.

The direct Coombs test demonstrates if the patient's RBCs have been attacked by antibodies in the patient's own bloodstream. Coombs serum is a solution containing antibodies to human globulin (antibodies). Coombs serum is mixed with the patient's RBCs. If the RBCs have antibodies on them, agglutination of the patient's RBCs will occur. The greater the quantity of antibodies against RBCs, the more clumping occurs. This test is read as *positive* with clumping on a scale of micropositive to +4. If the RBCs are not coated with autoantibodies against RBCs (immunoglobulins), agglutination will not occur; this is a *negative* test result.

### Interfering factors

- Antiphospholipid antibodies (see p. 65, anticardiolipin antibodies) can cause a false-positive DAT.
- ✚ Drugs that may cause false-positive results include ampicillin, captopril, cephalosporins, chlorpromazine, chlorpropamide, hydralazine, indomethacin, insulin, isoniazid, levodopa,

## 288 Coombs test, direct

methyldopa, penicillin, phenytoin, procainamide, quinidine, quinine, rifampin, streptomycin, sulfonamides, and tetracyclines.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or lavender
- Use venous blood from the umbilical cord to detect the presence of antibodies in the newborn
- List on the laboratory slip all medications and any transfusions that the patient has had in the last few days.

### Abnormal findings

Autoimmune hemolytic anemia

Transfusion reaction

Hemolytic disease of the newborn

Lymphoma

Systemic lupus erythematosus

Mycoplasma infection

Infectious mononucleosis

### notes

## Coombs test, indirect (Blood antibody screening, Indirect antiglobulin test [IAT])

**Type of test** Blood

### Normal findings

Negative; no agglutination

### Test explanation and related physiology

The indirect Coombs test detects circulating antibodies against RBCs. The major purpose of this test is to determine whether the patient has minor serum antibodies (other than the major ABO/Rh system) to RBCs that he or she is about to receive by blood transfusion. Therefore this test is the *screening* part of the *type and screen* routinely performed for blood compatibility testing (cross-matching in the blood bank).

Unlike the direct Coombs test that is performed on the patient's RBCs, this test is performed on the patient's serum. In this test, a small amount of the recipient's serum is added to donor RBCs containing known antigens on their surfaces. This is the first stage. In the second stage of the test, Coombs serum is added. If antibodies exist in the patient's serum, agglutination occurs. In blood transfusion screening, visible agglutination indicates that the recipient has antibodies to the donor's RBCs. If the recipient has no antibodies against the donor's RBCs, agglutination will not occur; transfusion should then proceed safely and without any transfusion reaction.

### Interfering factors

- ✚ Drugs that may cause false-positive results include antiarrhythmics, antituberculins, cephalosporins, chlorpromazine, insulin, levodopa, methyldopa, penicillin, phenytoin, quinidine, sulfonamides, and tetracyclines.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

Incompatible cross-matched blood

Maternal anti-Rh antibodies

Hemolytic disease of the newborn

Acquired immune hemolytic anemia

Presence of specific cold agglutinin antibody

### notes

**cortisol, blood, urine, saliva** (Hydrocortisone, Serum cortisol, Salivary cortisol, Total cortisol, Free cortisol)

**Type of test** Blood; urine; saliva

### Normal findings

#### Serum

Free cortisol

8 AM: 0.121-1.065 mcg/dL

Total cortisol:

Adult/elderly

8 AM: 5-23 mcg/dL or 138-635 nmol/L (SI units)

4 PM: 3-13 mcg/dL or 83-359 nmol/L (SI units)

Child 1-16 years

8 AM: 3-21 mcg/dL

4 PM: 3-10 mcg/dL

Newborn: 1-24 mcg/dL

#### Urine (24-hour)

Adult/elderly: < 100 mcg/24 hr or < 276 nmol/day (SI units)

Adolescent: 5-55 mcg/24 hr

Child: 2-27 mcg/24 hr

#### Saliva

7 AM-9 AM: 100-750 ng/dL

3 PM-5 PM: < 401 ng/dL

11 PM-midnight: < 100 ng/dL

### Test explanation and related physiology

The best method of evaluating adrenal activity is by directly measuring plasma cortisol levels. Normally, cortisol levels rise and fall during the day; this is called the diurnal variation. Cortisol levels are highest around 6 AM to 8 AM and gradually fall during the day, reaching their lowest point around midnight. Sometimes the earliest sign of adrenal hyperfunction is only the loss of this diurnal variation, even though the cortisol levels are not yet elevated. For example, individuals with Cushing syndrome often have upper normal plasma cortisol levels in the morning and do not exhibit a decline as the day proceeds. High levels of cortisol indicate Cushing syndrome, and low levels of plasma cortisol are suggestive of Addison disease.

For this test, blood is usually collected at 8 AM and again at around 4 PM. The 4 PM value is anticipated to be one-third to two-thirds of the 8 AM value. Normal values may be transposed in

individuals who have worked during the night and slept during the day for long periods of time.

The majority of cortisol circulates bound to corticosteroid-binding globulin (CBG) and albumin. Normally, less than 5% of circulating cortisol is free (unbound). Total cortisol includes measurements of free and bound cortisol.

The measurement of late-night *salivary cortisol* is another effective test for Cushing syndrome. It seems to be more convenient and superior to plasma and urine for detecting cortisol in patients with mild Cushing syndrome. Salivary cortisol assay cannot be used to diagnose hypocortisolism or Addison disease because laboratory methods are not sensitive enough at low levels. If late-night salivary cortisol levels are elevated, the results should be confirmed with a repeat salivary cortisol measurement, a midnight blood sampling for cortisol, or a 24-hour urinary collection of free cortisol. A dexamethasone suppression test (see p. 324) is another confirmation test that can be used.

### Interfering factors

- Pregnancy is associated with increased levels.
- Physical and emotional stress can elevate cortisol levels.
- Variations in protein levels caused by renal or liver disease can affect free cortisol levels.
- ✚ Drugs that may cause *increased* levels include amphetamines, cortisone, estrogen, oral contraceptives, and spironolactone (Aldactone).
- ✚ Drugs that may cause *decreased* levels include androgens, aminoglutethimide, betamethasone and other exogenous steroid medications, danazol, lithium, levodopa, metyrapone, and phenytoin (Dilantin).

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient to minimize anxiety.
- Assess the patient for signs of physical stress (e.g., infection) and report these to the physician.

#### During

##### *Blood*

- Collect a venous blood sample in a red-top or green-top tube in the morning after the patient has had a good night's sleep.
- Collect another blood sample at about 4 PM.
- Indicate the time of the venipuncture on the laboratory slip.

##### *Saliva*

1. Do not brush teeth before specimen collection.

## 292 cortisol, blood, urine, saliva

2. Do not eat or drink for 15 minutes before specimen collection.
3. Collect the specimen between 11 PM and midnight and record the collection time.
4. Collect at least 1.5 mL of saliva in a Salivette as follows:
  - a. Place swab directly into mouth by tipping container so that swab falls into mouth. Do not touch swab with fingers.
  - b. Keep swab in mouth for approximately 2 minutes. Roll swab in mouth; do not chew swab.
  - c. Place swab back into its container without touching and replace the cap.

### *Urine*

- PT** Instruct the patient how to collect a 24-hour urine. [See inside front cover for Routine Urine Testing.](#)
- Keep the collection on ice and use a preservative.

### **After**

- Apply pressure or a pressure dressing to the venipuncture site.

### **Abnormal findings**

#### ▲ **Increased levels**

Cushing syndrome  
Adrenal adenoma or carcinoma  
Ectopic ACTH-producing tumors  
Hyperthyroidism  
Obesity  
Stress

#### ▼ **Decreased levels**

Congenital adrenal hyperplasia  
Addison disease  
Hypopituitarism  
Hypothyroidism  
Liver disease

### **notes**



## C-peptide (Connecting peptide insulin, Insulin C-peptide, Proinsulin C-peptide)

**Type of test** Blood

### Normal findings

Fasting: 0.78-1.89 ng/mL or 0.26-0.62 nmol/L (SI units)

1 hour after glucose load: 5-12 ng/mL

### Test explanation and related physiology

In general, C-peptide levels correlate with insulin levels in the blood. The capacity of the pancreatic beta cells to secrete insulin can be evaluated by directly measuring either insulin or C-peptide. In most cases, direct measurement of insulin is more accurate. C-peptide levels, however, more accurately reflect islet cell function in the following situations:

- Patients with diabetes who are treated with exogenous insulin and who have antiinsulin antibodies.
- Patients who secretly administer insulin to themselves (factitious hypoglycemia). Insulin levels will be elevated. Direct insulin measurement in these patients tends to be high because the insulin measured is the self-administered exogenous insulin. But C-peptide levels in that same specimen will be low because exogenously administered insulin suppresses endogenous insulin (and C-peptide) production.
- Patients with diabetes who are taking insulin. This is done to see if the patient with diabetes is in remission and may not need exogenous insulin.
- Distinguishing type 1 from type 2 diabetes. This is particularly helpful in people who are newly diagnosed with diabetes. A person whose pancreas does not make any insulin (type 1 diabetes) has low levels of insulin and C-peptide. A person with type 2 diabetes has a normal or high level of C-peptide.

Furthermore, C-peptide is used in evaluating patients who are suspected of having an insulinoma. In patients with an autonomous secreting insulinoma, C-peptide levels are high. C-peptide can also be used to monitor treatment for insulinoma. A rise in C-peptide levels indicates a recurrence or progression of the insulinoma. Likewise, some clinicians use C-peptide testing as an indicator of the adequacy of therapeutic surgical pancreatectomy in patients with pancreatic tumors. C-peptide also can be used to diagnose insulin resistance syndrome.

## 294 C-peptide

### Interfering factors

- Because the majority of C-peptide is degraded in the kidney, renal failure can cause increased levels.
- ▼ Drugs that may cause *increased* levels of C-peptide include oral hypoglycemic agents (e.g., sulfonylureas).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Insulinoma  
Renal failure  
Pancreas transplant  
Type 2 diabetes mellitus

#### ▼ Decreased levels

Factitious hypoglycemia  
Radical pancreatectomy  
Type 1 diabetes mellitus

notes

## C-reactive protein test (CRP, High-sensitivity C-reactive protein [hs-CRP], Ultra-sensitive CRP)

**Type of test** Blood

### Normal findings

< 1.0 mg/dL or < 10.0 mg/L (SI units)

Cardiac risk:

Low: < 1.0 mg/dL

Average: 1.0-3.0 mg/dL

High: > 3.0 mg/dL

hs-CRP: < 3 mg/L

### Test explanation and related physiology

CRP is a nonspecific, acute-phase reactant used to diagnose bacterial infectious disease and inflammatory disorders, such as acute rheumatic fever and rheumatoid arthritis. CRP levels do not consistently rise with viral infections. CRP is a protein produced primarily by the liver during an acute inflammatory process and other diseases. A positive test result indicates the presence but not the cause of the disease. The synthesis of CRP is initiated by antigen-immune complexes, bacteria, fungi, and trauma.

The CRP test is a more sensitive and rapidly responding indicator than the erythrocyte sedimentation rate (ESR, see p. 373). In an acute inflammatory change, CRP shows an earlier and more intense increase than ESR; with recovery, the disappearance of CRP precedes the return of ESR to normal.

The level of CRP correlates with peak levels of the muscle/brain (MB) isoenzyme of creatine kinase (see p. 297), but CRP peaks occur 1 to 3 days later. Failure of CRP to normalize may indicate ongoing damage to the heart tissue. The level of CRP is a stronger predictor of cardiovascular events than the LDL cholesterol level. However, when used together with the lipid profile (in cholesterol, see p. 235), it adds prognostic information to that conveyed by the Framingham risk score.

The recent development of an assay for *high-sensitivity CRP* (*hs-CRP*) has enabled accurate assays at even low levels. Because of the individual variability in hs-CRP, two separate measurements are required to classify a person's risk level. In patients with stable coronary disease or acute coronary syndromes, hs-CRP measurement may be useful as an independent marker for assessing the likelihood of harmful events, including death, myocardial infarction, or restenosis after percutaneous coronary intervention.

### Interfering factors

- Elevated test results can occur in patients with hypertension, elevated body mass index, metabolic syndrome or diabetes mellitus, chronic infection (e.g., gingivitis, bronchitis), chronic inflammation (e.g., rheumatoid arthritis), and low HDL or high triglycerides.
- Cigarette smoking can cause increased levels.
- Decreased test levels can result from moderate alcohol consumption, weight loss, and increased activity/exercise.
- ✚ Drugs that may cause *increased* test results include estrogens and progesterones.
- ✚ Drugs that may cause *decreased* test results include fibrates, niacin, and statins.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: verify with laboratory
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

- Arthritis
- Acute rheumatic fever
- Reiter syndrome
- Crohn disease
- Vasculitis syndrome
- Systemic lupus erythematosus
- Tissue infarction or damage
- Acute myocardial infarction
- Pulmonary infarction
- Kidney transplant rejection
- Bone marrow transplant rejection
- Soft tissue trauma
- Bacterial infection
- Postoperative wound infection
- Urinary tract infection
- Tuberculosis
- Malignant disease
- Bacterial meningitis

**creatine kinase (CK, Creatine phosphokinase [CPK])****Type of test** Blood**Normal findings****Total CPK**

Adult/elderly

Male: 55-170 units/L or 55-170 units/L (SI units)

Female: 30-135 units/L or 30-135 units/L (SI units)

(Values are higher after exercise.)

Newborn: 68-580 units/L (SI units)

**Isoenzymes**

CK-MM: 100%

CK-MB: 0%


CK-BB: 0%

**Test explanation and related physiology**

This test is used to support the diagnosis of myocardial muscle injury (infarction). It can also indicate neurologic or skeletal muscle diseases. CK is found predominantly in the heart muscle, skeletal muscle, and brain. Serum CK levels are elevated when these muscle or nerve cells are injured. CK levels can rise within 6 hours after damage. If damage is not persistent, the levels peak at 18 hours after injury and return to normal in 2 to 3 days (Figure 12).

To test specifically for myocardial muscle injury, three CK isoenzymes are measured: CK-BB (CK1), CK-MB (CK2), and CK-MM (CK3). CK-MB appears to be specific for myocardial cells. CK-MB levels rise 3 to 6 hours after infarction occurs. If there is no further myocardial damage, the levels peak at 12 to 24 hours and return to normal 12 to 48 hours after infarction. CK-MB levels do not usually rise with transient chest pain caused by angina, pulmonary embolism, or congestive heart failure. One can expect to see a rise in CK-MB in patients with shock, malignant hyperthermia, myopathies, or myocarditis. Mild elevation of CK-MB (below the threshold of positive) can occur in patients with unstable angina and will signify an increased risk for an occlusive event. Very small amounts of CK-MB also exist in skeletal muscle. Severe injury to, or diseases of, the skeletal muscle can also raise the CK-MB above normal.

The CK-MB level is helpful in both quantifying the degree of myocardial infarction and timing the onset of infarction. The CK-MB level is often used to determine appropriateness

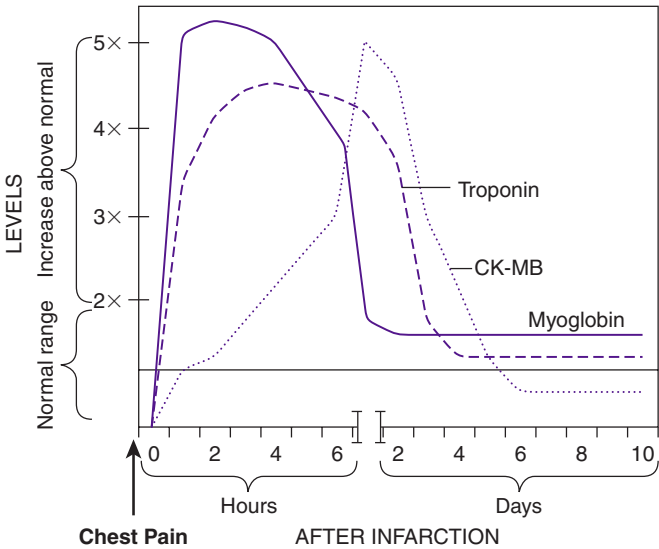

**FIG. 12** Blood studies useful in the diagnosis of myocardial infarction.

of thrombolytic therapy, which is used for myocardial infarction. High CK-MB levels suggest that significant infarction has already occurred, thereby precluding the benefit of thrombolytic therapy.

Because CK-BB is found predominantly in the brain and lung, injury to either of these organs (e.g., cerebrovascular accident, pulmonary infarction) will be associated with elevated levels of this isoenzyme.

CK-MM normally makes up almost all of the circulatory total CK enzymes in healthy people. When the total CK level is elevated as a result of increases in CK-MM, injury to or disease of the skeletal muscle is present. Examples of this include myopathies, vigorous exercise, multiple intramuscular (IM) injections, electroconvulsive therapy, cardioversion, chronic alcoholism, or surgery. Because CK is made only in the skeletal muscle, the normal value of total CK (and therefore CK-MM) varies according to a person's muscle mass. Large, muscular people may normally

have a CK level in the high range of normal. Likewise, people of small stature or those with low muscle mass will be expected to have low CK levels. This is important because high normal CK levels in these patients can mask an MI.

Each isoenzyme has been found to have isoforms. The CK-MM isoforms MM1 and MM3 are most useful for identifying cardiac disease. An MM3/MM1 ratio of greater than 1 suggests acute myocardial injury. A CK-MB ratio of MB2/MB1 greater than 1 also indicates acute myocardial injury.

CK is studied in patients with heart disease. Because its blood clearance and metabolism are well known, its frequent determination (on admission, at 12 hours, and at 24 hours) can accurately reflect timing, quantity, and resolution of an MI. Troponin (see p. 912) and myoglobin (see p. 635) are also serum markers used to confirm an MI (see [Figure 12](#)). A new assay is ischemia-modified albumin (see p. 550).

### Interfering factors

- IM injections may cause elevated CPK levels.
- Strenuous exercise and recent surgery may cause increased levels.
- Early pregnancy may cause decreased levels.
- ✚ Drugs that may cause *increased* levels include alcohol, amphotericin B, ampicillin, some anesthetics, anticoagulants, aspirin, captopril, colchicine, dexamethasone, fibrates, furosemide, lidocaine, lithium, morphine, propranolol, statins, and succinylcholine.

### Procedure and patient care

- See [inside front cover for Routine Blood Testing](#).
- Fasting: no
- Blood tube commonly used: red
- PT Discuss with the patient the need and reason for frequent venipuncture in diagnosing myocardial infarction.
  - Avoid IM injections in patients with cardiac disease. These injections may falsely elevate the total CPK level.
  - Blood collection is usually done daily for 3 days and then at 1 week.
  - Record the exact time and date of venipuncture on each laboratory slip. This aids in the interpretation of the temporal pattern of enzyme elevations.

## Abnormal findings

### ▲ Increased levels of total CK

Disease or injury affecting the heart muscle, skeletal muscle, or brain

### ▲ Increased levels of CK-BB isoenzyme

Disease affecting the central nervous system  
Adenocarcinoma (especially breast and lung)  
Pulmonary infarction

### ▲ Increased levels of CK-MB isoenzyme

Acute myocardial infarction  
Cardiac aneurysm surgery  
Cardiac defibrillation  
Myocarditis  
Ventricular arrhythmias  
Cardiac ischemia

### ▲ Increased levels of CK-MM isoenzyme

Rhabdomyolysis  
Muscular dystrophy  
Myositis  
Recent surgery  
Electromyography  
IM injections  
Crush injuries  
Delirium tremens  
Malignant hyperthermia  
Recent convulsions  
Electroconvulsive therapy  
Shock  
Hypokalemia  
Hypothyroidism  
Trauma

## notes



**creatinine, blood** (Serum creatinine)**Type of test** Blood**Normal findings**

Adult:

Female: 0.5-1.1 mg/dL or 44-97  $\mu\text{mol/L}$  (SI units)Male: 0.6-1.2 mg/dL or 53-106  $\mu\text{mol/L}$  (SI units)

Elderly: decrease in muscle mass may cause decreased values

Adolescent: 0.5-1.0 mg/dL

Child: 0.3-0.7 mg/dL

Infant: 0.2-0.4 mg/dL

Newborn: 0.3-1.2 mg/dL

**Possible critical values**

&gt; 4 mg/dL (indicates serious impairment in renal function)

**Test explanation and related physiology**

This test measures the amount of creatinine in the blood. Creatinine is a catabolic product of CPK, which is used in skeletal muscle contraction. The daily production of creatine, and subsequently creatinine, depends on muscle mass, which fluctuates very little. Creatinine, as with BUN (see p. 155), is excreted entirely by the kidneys and therefore is directly proportional to renal excretory function. Besides dehydration, only such renal disorders as glomerulonephritis, pyelonephritis, acute tubular necrosis, and urinary obstruction will cause abnormal elevations in creatinine.

The serum creatinine test, as with BUN, is used to diagnose impaired renal function. The creatinine test is used as an approximation of *glomerular filtration rate* (GFR). The serum creatinine level has much the same significance as the BUN level but tends to rise later. Therefore elevations in creatinine suggest chronicity of the disease process. In general, a doubling of creatinine suggests a 50% reduction in GFR. The creatinine level is interpreted in conjunction with the BUN test. These tests are referred to as *renal function studies*. The BUN/creatinine ratio is a good measurement of kidney and liver function. The normal adult range is 6 to 25, with 15.5 being the optimal adult value for this ratio.

Although serum creatinine is the most commonly used biochemical parameter to estimate GFR in routine practice, there are some shortcomings. Such factors as muscle mass and protein intake can influence serum creatinine, leading to an inaccurate

## 302 creatinine, blood

estimation of GFR. On the other hand, *cystatin C*, a protein that is produced at a constant rate by all nucleated cells, is probably a better indicator of GFR. Because of its constant rate of production, its serum concentration is determined only by glomerular filtration. Its level is not influenced by those factors that affect creatinine and BUN.

*Cystatin C* might predict the risk for developing *chronic kidney disease*, thereby signaling a state of preclinical kidney dysfunction. Several studies have found that increased levels of cystatin C are associated with the risk of death and several types of cardiovascular disease. For women, the average reference interval is 0.52 to 0.90 mg/L with a mean of 0.71 mg/L. For men, the average reference interval is 0.56 to 0.98 mg/L with a mean of 0.77 mg/L.

### Interfering factors

- ▼ Drugs that may *increase* creatinine values include ACE inhibitors, aminoglycosides, cimetidine, heavy-metal chemotherapeutic agents, and other nephrotoxic drugs (e.g., cephalosporins).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- For pediatric patients, blood is drawn from a heel stick.

### Abnormal findings

#### ▲ Increased levels

Glomerulonephritis  
Pyelonephritis  
Acute tubular necrosis  
Urinary tract obstruction  
Reduced renal blood flow  
(e.g., shock, dehydration,  
congestive heart failure,  
atherosclerosis)  
Diabetic nephropathy  
Nephritis  
Rhabdomyolysis  
Acromegaly  
Gigantism

#### ▼ Decreased levels

Debilitation  
Decreased muscle mass  
(e.g., muscular dystrophy,  
myasthenia gravis)

## creatinine clearance (CC, Estimated glomerular filtration rate [eGFR])

**Type of test** Urine (24-hour); blood

### Normal findings

Adult (< 40 years)

Male: 107-139 mL/min or 1.78-2.32 mL/s (SI units)

Female: 87-107 mL/min or 1.45-1.78 mL/s (SI units)

Newborn: 40-65 mL/min

Values decrease 6.5 mL/min/decade of life because of decline in GFR.

eGFR: > 60 mL/min/1.73 m<sup>2</sup>

### Test explanation and related physiology

Creatinine is a catabolic product of creatine phosphate, which is used in skeletal muscle contraction. The daily production of creatinine depends on muscle mass, which fluctuates very little. Creatinine is entirely excreted by the kidneys and therefore is directly proportional to the glomerular filtration rate (GFR; i.e., the number of milliliters filtered by all the nephrons in the kidneys per minute). The creatinine clearance (CC) is a measure of the GFR.

The CC depends on the amount of blood present to be filtered and the ability of the nephron to act as a filter. The amount of blood present for filtration is decreased in renal artery atherosclerosis, dehydration, or shock. The ability of the nephron to act as a filter is decreased by such diseases as glomerulonephritis, acute tubular necrosis, and most other primary renal diseases.

When one kidney becomes diseased, the opposite kidney, if normal, has the ability to compensate by increasing its filtration rate. Therefore with unilateral kidney disease or nephrectomy, a decrease in CC is not expected if the other kidney is normal.

Several nonrenal factors may influence CC. With each decade of age, the CC decreases 6.5 mL/min because of a decrease in the GFR. Urine collections are timed, and incomplete collections will falsely decrease CC. Muscle mass varies among people. Decreased muscle mass will give lower CC values.

The CC test requires a 24-hour urine collection and a serum creatinine level. The *uncorrected* CC is then computed using the following formula:

$$\text{Creatinine clearance} = \frac{UV}{P}$$

### 304 creatinine clearance

where

U = Number of milligrams per deciliter of creatinine excreted in the urine over 24 hours

V = Volume of urine in milliliters per minute

P = Serum creatinine in milligrams per deciliter

The *corrected* CC calculation takes into account the average body surface area.

The 24-hour urine collections used to measure CC are too time-consuming and expensive for routine clinical use. The GFR can be estimated (*eGFR*) using the Modification of Diet in Renal Disease (MDRD) Study equation. This is an equation that uses the serum creatinine, age, and numbers that vary depending on sex and ethnicity to calculate the GFR with very good accuracy. The prediction equation for GFR is as follows, with P<sub>cr</sub> being serum or plasma creatinine in mg/dL:

The GFR is expressed in mL/min/1.73 m<sup>2</sup>.

More and more, institutions across the country are beginning to report an eGFR on patients 18 years and older with every serum creatinine ordered. The eGFR calculation can be programmed into most laboratory information systems. As a result, chronic renal disease is being recognized more frequently in its early stages. Chronic kidney disease can be treated and progression to renal failure slowed or prevented. The eGFR can also be used to calculate medication dosage in patients with decreased renal function.

Table 13 shows population estimates for mean (average) eGFR by age. There is no difference between races or sexes when eGFRs are expressed per meter squared body surface area. For diagnostic purposes, most laboratories report eGFR values above 60 as greater than 60 mL/min/1.73 m<sup>2</sup>, not as an exact number.

**TABLE 13** Mean estimated glomerular filtration rates (eGFRs)

Age (years)	Mean eGFR (mL/min/1.73 m <sup>2</sup> )
20-29	116
30-39	107
40-49	99
50-59	93
60-69	85
70+	75

## Interfering factors

- Exercise may cause increased creatinine values.
- Incomplete urine collection may give a falsely lowered value.
- Pregnancy increases CC.
- A diet high in meat can transiently elevate CC.
- The eGFR may be inaccurate in extremes of age and in patients with obesity, severe malnutrition, paraplegia, quadriplegia, or pregnancy.
- ▼ Drugs that may cause *increased* levels include aminoglycosides, cimetidine, heavy-metal chemotherapeutic agents, and nephrotoxic drugs such as cephalosporins.
- ▼ Drugs that may cause *decreased* eGFR are drugs that interfere with creatinine secretion (e.g., cimetidine or trimethoprim) or creatinine assay (cephalosporins).

## Procedure and patient care

- See inside front cover for Routine Urine Testing.
- PT Note that some laboratories instruct the patient to avoid cooked meat, tea, coffee, or drugs on the day of the test.
- Make sure a venous blood sample is drawn in a red-top tube during the 24-hour urine collection.
- Mark the patient's age, weight, and height on the requisition sheet.

## Abnormal findings

### ▲ Increased levels

Exercise  
Pregnancy  
High cardiac output syndromes

### ▼ Decreased levels

Impaired kidney function  
(e.g., renal artery atherosclerosis, glomerulonephritis, acute tubular necrosis)  
Conditions causing decreased GFR (e.g., congestive heart failure, cirrhosis with ascites, shock, dehydration)

notes

## cryoglobulin

**Type of test** Blood

### Normal findings

No cryoglobulins detected

### Test explanation and related physiology

Cryoglobulins are abnormal globulin protein complexes that exist in the blood of patients with various diseases. These proteins precipitate reversibly at low temperatures and redissolve with rewarming. They can precipitate in the blood vessels of the fingers when exposed to cold temperatures. This precipitation causes sludging of the blood in those blood vessels. These patients may have symptoms of purpura, arthralgia, or Raynaud phenomenon.

Type I (monoclonal) cryoglobulinemia is associated with monoclonal gammopathy of undetermined significance, macroglobulinemia, or multiple myeloma. Type II (mixed, two or more immunoglobulins of which one is monoclonal) cryoglobulinemia is associated with autoimmune disorders, such as systemic lupus erythematosus, rheumatoid arthritis, and Sjögren syndrome. Type III (polyclonal) cryoglobulinemia is associated with the same disease spectrum as type II cryoglobulinemia.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: verify with laboratory
- Blood tube commonly used: red

**PT** Inform the patient that an 8-hour fast may be required.

**PT** If cryoglobulins are found to be present, warn the patient to avoid cold temperatures and contact with cold objects to minimize Raynaud symptoms.

### Abnormal findings

Connective tissue disease (e.g., lupus erythematosus, Sjögren syndrome, rheumatoid arthritis)

Lymphoid malignancies (e.g., multiple myeloma, leukemia, Waldenström macroglobulinemia, lymphoma)

Acute and chronic infections (e.g., infectious mononucleosis, endocarditis, poststreptococcal glomerulonephritis)

Liver disease (e.g., hepatitis, cirrhosis)

### notes

## **cystography** (Cystourethrography, Voiding cystography, Voiding cystourethrography [VCUG])

**Type of test** X-ray with contrast dye

### **Normal findings**

Normal bladder structure and function

### **Test explanation and related physiology**

Filling the bladder with radiopaque contrast material provides visualization of the bladder for radiographic study. Either fluoroscopic or x-ray images demonstrate bladder filling and collapse after emptying. Filling defects or shadows within the bladder indicate primary bladder tumors. Extrinsic compression or distortion of the bladder is seen with pelvic tumor (e.g., rectal, cervical) or hematoma (secondary to pelvic bone fractures). Extravasation of the dye is seen with traumatic rupture, perforation, and fistula of the bladder. Vesicoureteral reflux (abnormal backflow of urine from bladder to ureters), which can cause persistent or recurrent pyelonephritis, also may be demonstrated during cystography. Although the bladder is visualized during an IV pyelogram (see p. 764), primary pathologic bladder conditions are best studied by cystography.

### **Contraindications**

- Patients with urethral or bladder infection or injury

### **Potential complications**

- Urinary tract infection
- This may result from catheter placement or the instillation of contaminated contrast material.
- Allergic reaction to iodinated dye. This rarely occurs because the dye is not administered intravenously.

### **Procedure and patient care**

#### **Before**

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- Obtain informed consent if required by the institution.
  - Give clear liquids for breakfast on the morning of the test.
- PT** Assure the patient that he or she will be draped to prevent unnecessary exposure.
- Insert a Foley catheter if ordered.

## 308 cystography

### During

- Note the following procedural steps:
    1. The patient is taken to the radiology department and placed in a supine or lithotomy position.
    2. Unless the catheter is already present, one is placed.
    3. Through the catheter, approximately 300 mL of air or radiopaque dye (much less for children) is injected into the bladder.
    4. The catheter is clamped.
    5. X-ray images are taken.
    6. If the patient is able to void, the catheter is removed, and the patient is asked to urinate while images are taken of the bladder and urethra (voiding cystourethrogram). Ensure that males wear a lead shield over the testes to prevent irradiation of the gonads.
  - Remember that female patients cannot be shielded without blocking bladder visualization.
  - Note that a radiologist performs the study in approximately 15 to 30 minutes.
- PT** Tell the patient that this test is moderately uncomfortable if bladder catheterization is required.

### After

- Assess the patient for signs of urinary tract infection.
- PT** Encourage the patient to drink fluids to eliminate the dye and to prevent accumulation of bacteria.

### Abnormal findings

Bladder tumor

Pelvic tumor

Hematoma

Bladder trauma

Vesicoureteral reflux

### notes




**cystoscopy** (Endourology)**Type of test** Endoscopy**Normal findings**

Normal structure and function of the urethra, bladder, ureters, and prostate (in males)

**Test explanation and related physiology**

Cystoscopy provides direct visualization of the urethra and bladder through the transurethral insertion of a cystoscope into the bladder (Figure 13). Cystoscopy is used *diagnostically* to allow

- Direct inspection and biopsy of the prostate, bladder, and urethra
- Collection of a separate urine specimen directly from each kidney by the placement of ureteral catheters

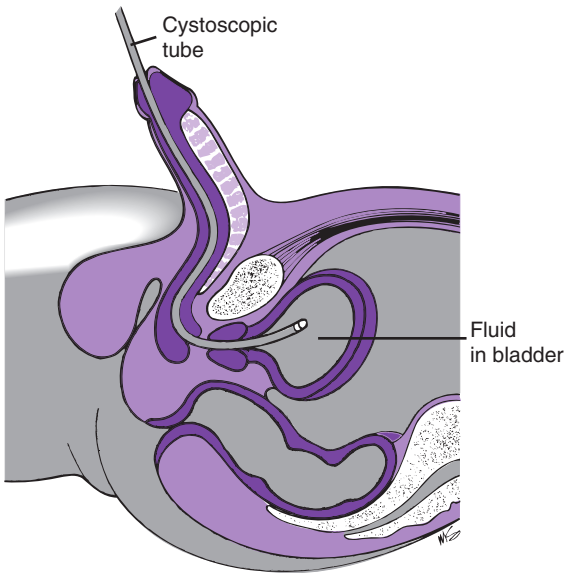


**FIG. 13** Cystoscopic examination of the male bladder. The cystoscope is passed through the urethra into the bladder. Although shown here as a flexible scope, usually the scope is rigid. Through the scope, fluid is instilled to maintain bladder distention.

### 310 cystoscopy

- Measurement of bladder capacity and determination of ureteral reflux
- Identification of bladder and ureteral calculi
- Placement of ureteral catheters (Figure 14) for retrograde pyelography (see p. 310)
- Identification of the source of hematuria

Cystoscopy is used *therapeutically* to provide

- Resection of small, superficial bladder tumors
- Removal of foreign bodies and stones
- Dilatation of the urethra and ureters
- Placement of stents to drain urine from the renal pelvis
- Coagulation of bleeding areas
- Implantation of radium seeds into a tumor
- Resection of hypertrophied or malignant prostate gland overgrowth
- Placement of ureteral stents during pelvic surgery

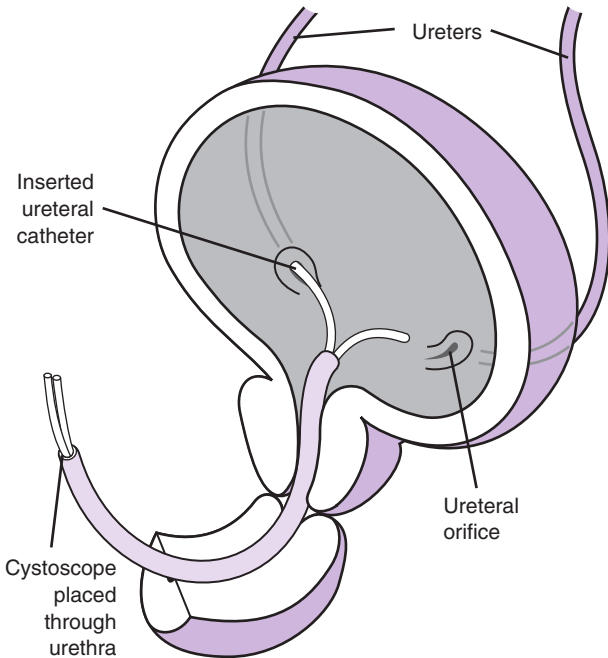

**FIG. 14** Ureteral catheterization through the cystoscope. Note the ureteral catheter inserted into the right orifice. The left ureteral catheter is ready to be inserted.

The cystoscope consists primarily of an obturator and a telescope. The obturator is used to insert the cystoscope atraumatically. After the cystoscope is in the bladder, the obturator is removed, and the telescope is passed through the cystoscope. The lens and lighting system of the telescope permit adequate visualization of the lower genitourinary tract. Transendoscopic instruments (e.g., forceps, scissors, needles, electrodes) are used when appropriate. *Endourology* refers to endoscopic surgery performed on the bladder and urethra during cystoscopy.

### Potential complications

- Perforation of the bladder
- Sepsis caused by bacteria from infected urine
- Hematuria
- Urinary retention

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Ensure that an informed consent is obtained.
- If enemas are ordered to clear the bowel, assist the patient as needed.

**PT** Encourage the patient to drink fluids several hours before the procedure to maintain a continuous flow of urine for collection and to prevent multiplication of bacteria.

- If the procedure will be done with the patient under local anesthesia, allow a liquid breakfast.
- If the procedure will be performed under general anesthesia, follow routine precautions. Keep the patient NPO after midnight on the day of the test.
- Administer the preprocedure medications as ordered before the study. Sedatives decrease the spasm of the bladder sphincter, decreasing the patient's discomfort.

#### During

- Note the following procedural steps:
  1. Cystoscopy is performed in the operating room or in the urologist's office.
  2. The patient is placed in the lithotomy position with his or her feet in stirrups.
  3. The external genitalia are cleansed with an antiseptic.
  4. A local anesthetic is instilled into the urethra if general anesthesia has not been used.
  5. The cystoscope is inserted, and the desired diagnostic or therapeutic studies are performed.

## 312 cystoscopy

- PT** Instruct the patient to lie very still during the entire procedure to prevent trauma to the urinary tract.
- PT** Tell the patient that he or she will have the desire to void as the cystoscope passes the bladder neck.
  - When the procedure is completed, keep the patient on bed rest for a short time.
  - Note that if endourology is performed, the urethra will also be evaluated.
  - Note that this procedure is performed by a urologist in approximately 25 minutes.
- PT** When local anesthesia is used, inform the patient of the associated discomfort (much more than with urethral catheterization).

### After

- PT** Instruct the patient not to walk or stand alone immediately after his or her legs have been removed from the stirrups. Orthostasis may result from standing erect.
  - Assess the patient's ability to void for at least 24 hours after the procedure. Urinary retention may be secondary to edema caused by instrumentation.
- PT** Instruct the patient to note the urine color. Pink-tinged urine is common. The presence of bright red blood or clots should be reported to the physician.
  - Monitor the patient for complaints of back pain, bladder spasms, urinary frequency, and burning on urination. Warm sitz baths and mild analgesics may be ordered and given. Sometimes belladonna and opium (B&O) suppositories are given to relieve bladder spasms. Warm, moist heat to the lower abdomen may help relieve pain and promote muscle relaxation.
- PT** Encourage increased intake of fluids. A dilute urine decreases dysuria. Fluids also maintain a constant flow of urine to prevent stasis and the accumulation of bacteria in the bladder.
  - Monitor the patient's vital signs for a decrease in blood pressure and an increase in pulse as an indication of hemorrhage.
  - Observe for signs and symptoms of sepsis (elevated temperature, flush, chills, decreased blood pressure, increased pulse).
- PT** Note that antibiotics are occasionally ordered before and after the procedure to reduce the incidence of bacteremia that may result from instrumentation.
- PT** Instruct the patient to watch for fever, shaking chills, or prolonged dysuria as possible signs of a urinary tract infection.

- PT** Encourage the patient to use cathartics, especially after cystoscopic surgery. Increases in intraabdominal pressure caused by constipation may initiate urologic bleeding.
- If postprocedure irrigation is ordered, use an isotonic solution containing mannitol, glycine, or sorbitol to prevent fluid overhydration in the event any of the irrigation is absorbed through opened venous sinuses in the bladder.
- PT** If a catheter is left in after the procedure, provide catheter care instructions.

### Abnormal findings

Lower urologic tract tumor  
Stones in the ureter or bladder  
Prostatic hypertrophy  
Prostate cancer  
Inflammation of the bladder and urethra  
Urethral or ureteral stricture  
Prostatitis  
Vesical neck contracture

### notes

**cytokines** (interleukins)**Type of test** Blood**Normal findings**

Vary by laboratory and technique

**Test explanation and related physiology**

Cytokines were renamed as *interleukins* and were numbered by the sequence of discovery. Interleukins, in general, are made up of leukocytes. Lymphokines and monokines are made up of lymphocytes and monocytes, respectively. Other cytokines include *interferon* and *growth factors*.

Some cytokines are produced at increased levels in certain disease states and therefore become markers for disease extent, progression, and response to therapy. They can act as tumor markers in cancers associated with elevated cytokines. *Human interferon inducible protein 10* is a small cytokine belonging to the chemokine family that affects cellular chemotaxis, immune response, and bone marrow inhibition. This protein, when present in high quantities in an acutely ill patient, is an accurate predictor of multiple organ failure.

Any table created to list all the cytokines and their functions quickly becomes inaccurate. The discovery of new cytokines and new functions changes so frequently that any such table would be outdated in the delay to publication. Likewise, any listing of normal values would be just as quickly obsolete because methods of testing change so frequently. It is suggested that reference to normal values be directed to the laboratory performing the assay.

At present, cytokine quantitative and qualitative assays are predominantly used for research. Clinically, cytokine assays may have the following uses:

- Measurement of AIDS progression
- Measurement of progression of inflammatory diseases, such as rheumatoid arthritis
- Tumor markers (e.g., breast cancer, lymphoma, and leukemia)
- Determination of disease risk (e.g., risk of developing Kaposi sarcoma in AIDS patients)
- Determination of treatment of disease (e.g., patients with rheumatoid arthritis)
- Determination of immune function and response
- Monitoring of patients receiving cytokine therapy or anticytokine therapy

Usually, cytokine testing is performed on serum. However, joint fluid is often tested in the evaluation of the patient with arthritis. Likewise, if inflammatory encephalitis or meningitis is considered, cerebrospinal fluid may be the specimen.

### Interfering factors

- Cells can still produce cytokines after specimen collection. It is best to freeze the specimen.
- Cytokines can degrade in the specimen container.
- Cytokines can stimulate or inhibit other cytokines while in the specimen container.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

Inflammatory disease

AIDS

Various malignancies

notes

**cytolethal distending toxin B (CdtB) and antivinculin antibodies****Type of test** Blood**Normal findings**

1-2.5 titers

**Test explanation and related physiology**

With no clear pathophysiology, irritable bowel syndrome (IBS) is currently identified through a diagnosis of exclusion and is considered a functional disease. The Rome IV criteria are used to diagnose IBS. These criteria include the following: symptoms of recurrent abdominal pain or discomfort with a marked change in bowel habit for at least 6 months, symptoms on at least 3 days for at least 3 months, and two or more of the following pain descriptions:

- Pain relieved by a bowel movement
- Pain onset linked to a change in frequency of stool
- Pain onset linked to a change in appearance of stool

Patients undergo extensive testing to rule out inflammatory bowel disease (IBD), and maldigestion (such as celiac disease). Diarrhea-predominant irritable bowel syndrome (D-IBS) can be precipitated by acute gastroenteritis. Anti-CdtB and antivinculin antibodies may occur after an acute bout of gastroenteritis. Persistent elevation of these antibodies is predictive of D-IBS rather than IBD or celiac disease.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: A central laboratory testing kit
- Obtain a list of foods that have been ingested in the last 48 hours.
- Assess severity of symptoms and the hydration status of the patient with diarrhea.

**Abnormal findings**

Diarrhea-predominate irritable bowel syndrome

Acute gastroenteritis

**notes**



## cytomegalovirus (CMV)

**Type of test** Blood

### Normal findings

No virus isolated

### Test explanation and related physiology

CMV is part of the viral family that includes herpes simplex, Epstein–Barr, and varicella-zoster viruses. CMV infection is widespread and common. Infections usually occur in the fetus, during early childhood, and in the young adult. Certain populations are at increased risk. Male homosexuals, transplant patients, and AIDS patients are particularly susceptible. Infections are acquired by contact with body secretions or urine. Blood transfusions are a common form of spread for CMV.

CMV is the most common congenital infection. Pregnant mothers can get the disease during pregnancy, or a previous CMV infection can become reactivated. Approximately 10% of infected newborns exhibit permanent damage, usually mental retardation and auditory damage. Fetal infection can cause microcephaly, hydrocephaly, cerebral palsy, mental retardation, or death. The term TORCH (*toxoplasmosis, other, rubella, cytomegalovirus, herpes*) has been applied to infections with recognized detrimental effects on fetuses.

Virus culture is the most definitive method of diagnosis. Antibodies reveal much more information about the activity of the infection. CMV IgG antibody levels persist for years after infection. Identification of IgM antibodies, however, indicates a relatively recent primary infection. Three different CMV antigens can be detected immunologically. They are called *early*, *intermediate-early*, and *late* antigens and indicate onset of infection. A fourfold increase in CMV titer in paired sera drawn 10 to 14 days apart is usually indicative of an acute infection.

More recently, measurement of *CMV-specific IgG avidity* is able to distinguish primary from nonprimary CMV infections. In this test, the strength with which the IgG attaches to the CMV antigen is measured. IgG avidity matures with the length of time after primary infection. Therefore IgG produced in the first few months after primary CMV infection will exhibit “low avidity.” IgG produced more than 6 to 8 months after CMV infection will have “high avidity” and represent nonprimary chronic CMV infection.

## 318 cytomegalovirus

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or gold
- For culture specimens, a urine, sputum, or mouth swab is the specimen of choice. Fresh specimens are essential.
- The specimens are cultured in a virus laboratory, which takes about 3 to 7 days.
- For antibody or antigen titer, collect a venous blood sample in a gold- or red-top tube.
- Collect a specimen from a mother with suspected acute infection as early as possible.
- Collect the convalescent specimen 2 to 4 weeks later.

### Abnormal findings

CMV infection

notes

**D-dimer** (Fragment D-dimer, Fibrin degradation product [FDP], Fibrin split products)**Type of test** Blood**Normal findings**

&lt; 250 ng/mL or &lt; 0.4 mcg/mL

**Test explanation and related physiology**

The fragment D-dimer test assesses both thrombin and plasmin activity. D-dimer is a fibrin degradation fragment that is made through lysis of crosslinked (D-dimerized) fibrin. As plasmin acts on the fibrin polymer clot, fibrin degradation products (FDPs) and D-dimer are produced. The D-dimer assay provides a highly specific measurement of the amount of fibrin degradation that occurs. Normal plasma does not have detectable amounts of fragment D-dimer. For discussion of other fibrin degradation products, see thrombosis indicators (see p. 877).

This test provides a simple and confirmatory test for disseminated intravascular coagulation (DIC). Positive results of the D-dimer assay correlate with positive results of other thrombosis indicators. The D-dimer assay may be more specific than the FDP assay, but it is less sensitive. Therefore combining the FDP and D-dimer assays provides a highly sensitive and specific test for recognizing DIC.

Levels of D-dimer can also increase when a fibrin clot is lysed by thrombolytic therapy. Thrombotic problems, such as deep vein thrombosis (DVT), pulmonary embolism (PE), sickle cell anemia, and thrombosis of malignancy, are associated with high D-dimer levels. D-dimer is used as an effective screening test for DVT. It is able to accurately identify patients with DVT who are then sent for venous duplex scanning (see p. 960). The D-dimer test, however, is often positive in patients who are already hospitalized. If the D-dimer test result is negative, its high predictability indicates that the patient does not have PE/DVT; further testing may not be necessary. D-dimer is an extremely sensitive and specific test for PE. Among patients presenting to the emergency department, normal D-dimer levels have a high negative predictive value for PE.

Finally, the D-dimer test can be used to determine the duration of anticoagulation therapy in patients with DVT. Patients with an abnormal D-dimer level 1 month after the discontinuation of anticoagulant therapy have a significant incidence of

## 320 D-dimer

recurrent DVT. This incidence can be reduced by restarting anti-coagulation therapy.

### Interfering factors

- D-dimer level may be decreased in lipemic patients.
- The presence of rheumatoid factor at a level  $> 50$  IU/mL may lead to increased D-dimer levels.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: blue
- Remember that if the patient is receiving anticoagulants or has coagulopathies, the bleeding time will be increased.

### Abnormal findings

#### ▲ Increased levels

Fibrinolysis

During thrombolytic or defibrination therapy with tissue plasminogen activator

Deep vein thrombosis

Pulmonary embolism

Arterial thromboembolism

Disseminated intravascular coagulation

Sickle cell anemia with or without vasoocclusive crisis

Pregnancy

Malignancy

Surgery

notes

**delta-aminolevulinic acid** (Aminolevulinic acid [ALA], d-ALA)**Type of test** Urine (24-hour)**Normal findings**1.5-7.5 mg/24hr or 11-57  $\mu\text{mol}/24\text{ hr}$  (SI units)**Possible critical values**

&gt; 20 mg/24hr

**Test explanation and related physiology**

As the basic precursor for the porphyrins (see p. 717), d-ALA is needed for the normal production of porphobilinogen, which ultimately leads to heme synthesis in erythroid cells. Heme is used in the synthesis of hemoglobin. Genetic disorders (porphyria) are associated with a lack of a particular enzyme vital to heme metabolism. These disorders are characterized by an accumulation of porphyrin products in the liver or red blood cells. Acute intermittent porphyria (AIP) is the most common form of the liver porphyrias.

In lead intoxication, heme synthesis is similarly diminished by the inhibition of ALA dehydrase. This enzyme assists in the conversion of ALA to porphobilinogen. As a result of lead poisoning, ALA accumulates in the blood and urine.

**Interfering factors**

✘ Drugs that may cause *increased* ALA levels include barbiturates, griseofulvin, and penicillin.

**Procedure and patient care**

- See inside front cover for Routine Urine Testing.
- Keep the urine in a light-resistant container with a preservative.
- If the patient has a Foley catheter in place, cover the drainage bag to prevent exposure to light.

**Abnormal findings****▲ Increased levels**

Porphyria  
Lead intoxication  
Chronic alcoholic disorders  
Diabetic ketoacidosis

**notes**

**dental x-ray** (Dental radiography)**Type of test** X-ray**Normal findings**

No evidence of dental caries (tooth decay/cavity), tumor, tooth impaction, or congenital abnormalities

**Test explanation and related physiology**

Dental x-rays can be taken intraorally and extraorally. Intraoral x-rays (with the film or digital plate inside the mouth) are the most common type of radiograph taken in dentistry. They give a high level of detail of the tooth, bone, and supporting tissues of the mouth. These x-rays allow dentists to:

- Find cavities.
- Look at the tooth roots.
- Check the health of the bony area around the tooth.
- Determine periodontal disease.
- See the status of developing teeth.

Intraoral x-rays show different aspects of the teeth. They include several different views:

- The bitewing view—highlights the crowns of the molars and bicuspids. These views find decay particularly between back teeth.
- The periapical view—highlights only one or two teeth at a time and images the entire length of each tooth, from crown to root.
- The occlusal view—highlights tooth development and placement in children. This view images nearly the full arch of teeth in either the upper or lower jaw.

Extraoral x-rays are made with the film or digital plate outside the mouth. They demonstrate the teeth, but they also provide information on the jaw and skull. They are used to:

- Keep track of growth and development.
- Look at the status of impacted teeth.
- Examine the relationships between teeth and jaws.
- Examine the bones of the face.

Extraoral images include panoramic x-rays, tomograms, cephalometric views, and CT scans. Extraoral x-rays are less detailed than intraoral x-rays. For this reason, they are not usually used for detecting cavities or flaws in individual teeth.

In most dental offices, images are still produced on x-ray film. However, many dental x-ray facilities are switching over to

digital radiography in which the images are viewed on a computer monitor.

### Interfering factors

- Earrings can obstruct radiographic visualization of part of the bone to be evaluated.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See page xviii for a discussion of radiation exposure and risk.
- PT** Instruct the patient that he or she will need to keep still while the x-ray image is being taken (about 1 second).
- Shield the patient's testes, ovaries, or pregnant abdomen to avoid exposure from scattered radiation.
- PT** Tell the patient that no fasting or sedation is required.

#### During

- Note that this test is routinely performed by a dental hygienist within several minutes.
- PT** Tell the patient that no discomfort is associated with this test, except possibly the large appliance for intraoral images.

#### After

- PT** Explain the radiographic interpretation to the patient.

### Abnormal findings

Dental caries

Dental malformation and growth

Bone or jaw tumors, resorption, infection

Tooth injury

Periodontal or endodontal disease

#### notes

**dexamethasone suppression test** (DST, Prolonged/rapid DST, Cortisol suppression test, ACTH suppression test)

**Type of test** Blood; urine (24-hour)

**Normal findings****Prolonged method**

Low dose: > 50% reduction of plasma cortisol

High dose: > 50% reduction of plasma cortisol

Urinary free cortisol: < 20 mcg/24 hr (< 50 nmol/24 hr)

**Rapid (overnight) method**

Normal: plasma cortisol levels suppressed to < 2 mcg/dL

**Test explanation and related physiology**

The DST is based on pituitary adrenocorticotropic hormone (ACTH) secretion being dependent on the plasma cortisol feedback mechanism. As plasma cortisol levels increase, ACTH secretion is suppressed; as cortisol levels decrease, ACTH secretion is stimulated. Dexamethasone is a synthetic steroid (similar to cortisol) that normally should suppress ACTH secretion. Under normal circumstances, this results in reduced stimulation to the adrenal glands and ultimately a drop of 50% or more in plasma cortisol and 17-OHCS levels. This important feedback system does not function properly in patients with hypercortisol states.

In Cushing syndrome caused by bilateral adrenal hyperplasia (Cushing disease), the pituitary gland is reset upward and responds only to high plasma levels of cortisol or its analogs. In Cushing syndrome caused by adrenal adenoma or cancer (which acts autonomously), cortisol secretion continues despite a decrease in ACTH. When Cushing syndrome is caused by an ectopic ACTH-producing tumor (as in lung cancer), that tumor is also considered autonomous and will continue to secrete ACTH despite high cortisol levels. Again, no decrease occurs in plasma cortisol. Knowledge of the following defects in the normal cortisol-ACTH feedback system is the basis for differentiating hypercortisol states using DST. ACTH and plasma cortisol levels are measured during this test.

When hypercortisol is caused by:

- *Bilateral adrenal hyperplasia (Cushing disease)*
  - Low dose: no change
  - High dose: > 50% reduction of plasma cortisol and ACTH is elevated



- *Adrenal adenoma or carcinoma (primary hypercortisolism)*
  - Low dose: no change
  - High dose: no change
  - ACTH is undetectable or low
- *Ectopic ACTH-producing tumor*
  - Low dose: no change
  - High dose: no change
  - ACTH is normal to elevated

The DST may also identify depressed persons likely to respond to electroconvulsive therapy or antidepressants rather than to psychological or social interventions. ACTH production will not be fully suppressed after administration of low-dose dexamethasone in these patients.

The *prolonged* DST can be performed over a 2-day period on an outpatient basis. The *rapid* DST is easily and quickly performed and is used primarily as a screening test to diagnose Cushing syndrome. It is less accurate and informative than the prolonged DST, but when its results are normal, the diagnosis of Cushing syndrome can be safely excluded.

### Interfering factors

- Physical and emotional stress can elevate ACTH release.
- ✦ Drugs that can affect test results include barbiturates, estrogens, oral contraceptives, phenytoin, spironolactone, steroids, and tetracyclines.

### Procedure and patient care

#### Before

**PT** Explain the procedure (prolonged or rapid test) to the patient.

#### During

- There are several documented methods of performing this test by varying the dose and duration of testing.

#### *Prolonged test*

- Obtain a baseline 24-hour urine collection for urinary free cortisol (see p. 290). *See inside front cover for Routine Urine Testing.*
- Collect blood for determination of baseline plasma cortisol levels (see p. 290) if indicated.
- Collect 24-hour urine specimens daily over a 2-day period. Because 2 continuous days of urine collections are needed, no urine specimens are discarded except for the first voided specimen on day 1, after which the collection begins.
- On days 1 and 2, administer a low dose of dexamethasone by mouth every 6 hours for 48 hours.

### 326 dexamethasone suppression test

- Administer the dexamethasone with milk or an antacid to prevent gastric irritation.
- Note that the urine samples for free cortisol do not need a preservative.
- Note that creatinine is measured in all the 24-hour urine collections to demonstrate their accuracy and adequacy.
- Keep the urine specimens refrigerated or on ice during the collection period.

#### *Rapid test*

- Give the patient a low dose of dexamethasone by mouth at 11 PM.
- Administer the dexamethasone with milk or an antacid.
- Attempt to ensure a good night's sleep. However, use sedative-hypnotic drugs only if absolutely necessary.
- At 8 AM the next morning, draw blood for determination of plasma cortisol level before the patient arises.
- If no cortisol suppression occurs after the dose of dexamethasone, administer a higher dose at 11 PM and obtain a cortisol level as described previously. This is referred to as the *overnight dexamethasone suppression test*. Patients with adrenal hyperplasia will suppress. Patients with adrenal or ectopic tumors will not suppress.

#### **After**

- Assess the patient for steroid-induced side effects by monitoring glucose levels and potassium levels.

#### **Abnormal findings**

Cushing syndrome

Cushing disease

Ectopic ACTH-producing tumors

Adrenal adenoma or carcinoma

Bilateral adrenal hyperplasia

Mental depression

Hyperthyroidism

#### **notes**

## diabetes mellitus autoantibody panel (Insulin autoantibody [IAA], Islet cell antibody [ICA], Glutamic acid decarboxylase antibody [GAD Ab])

**Type of test** Blood

### Normal findings

< 1:4 titer; no antibody detected

### Test explanation and related physiology

Type 1 diabetes mellitus (DM) is insulin-dependent (IDDM). It is becoming increasingly recognized that this disease is an *organ-specific* form of autoimmune disease that results in destruction of the pancreatic islet cells and their products. These antibodies are used to differentiate type I DM from type II non-insulin-dependent DM. Nearly 90% of people with type 1 DM have one or more of these autoantibodies at the time of their diagnosis. People with type 2 DM have low or negative titers.

These antibodies often appear years before the onset of symptoms. The panel is useful to screen relatives of IDDM patients who are at risk of developing the disease. GAD Ab provide confirmatory evidence. The presence of these antibodies identifies which woman with gestational diabetes will eventually require insulin permanently.

The presence of insulin antibodies is diagnostic of factitious hypoglycemia from surreptitious administration of insulin. This antibody panel is also used in surveillance of patients who have received pancreatic islet cell transplantation. Finally, these antibodies can be used to identify late onset type 1 diabetes in those patients previously thought to have type 2 diabetes.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or serum separator

### Abnormal findings

#### ▲ Increased levels

- Insulin-dependent diabetes mellitus
- Insulin resistance
- Allergies to insulin
- Factitious hypoglycemia

**2,3-diphosphoglycerate (2,3-DPG in erythrocytes)**

**Type of test** Blood

**Normal findings**

$12.3 \pm 1.87 \mu\text{mol/g}$  of hemoglobin or  $0.79 \pm 0.12 \text{ mol/mol}$  hemoglobin (SI units)

$4.2 \pm 0.64 \mu\text{mol/mL}$  of erythrocytes or  $4.2 \pm 0.64 \text{ mmol/L}$  erythrocytes (SI units)

(Levels are lower in newborns and even lower in premature infants.)

**Test explanation and related physiology**

This test is used in the evaluation of nonspherocytic hemolytic anemia. 2,3-DPG is a byproduct of the glycolytic respiratory pathway of the red blood cell (RBC). A congenital enzyme deficiency in this vital pathway alters the RBC shape and survival significantly. Anemia is the result. Another result of the enzyme deficiency is reduced synthesis of 2,3-DPG. Because 2,3-DPG controls oxygen transport from RBCs to tissues, deficiencies of this enzyme result in alterations of the RBC-oxygen dissociation curve that controls the release of oxygen to tissues. Many anemias not caused by 2,3-DPG deficiency are associated with increased levels of 2,3-DPG as a compensatory mechanism.

Usually, 2,3-DPG levels increase in response to anemia or hypoxic conditions (e.g., obstructive lung disease, congenital cyanotic heart disease, after vigorous exercise). Increases in 2,3-DPG decrease the oxygen binding to hemoglobin so that oxygen is more easily released to the tissues when needed (lower arterial  $\text{PO}_2$ ). Levels of 2,3-DPG are decreased as a result of inherited genetic defects. This genetic defect parallels that of sickle cell anemia and hemoglobin C diseases.

**Interfering factors**

- Vigorous exercise may cause increased levels.
- High altitudes may increase levels.
- Banked blood has decreased amounts of 2,3-DPG.
- Acidosis decreases levels.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

## Abnormal findings

### ▲ Increased levels

Anemia  
Hypoxic heart (e.g., cyanotic heart disease) or lung (e.g., chronic obstructive pulmonary disease) diseases  
Hyperthyroidism  
Chronic renal failure  
Pyruvate kinase deficiency  
Cystic fibrosis  
Adjustment to higher altitudes

### ▼ Decreased levels

Polycythemia  
Acidosis  
Post-massive blood transfusion  
2,3-DPG mutase deficiency  
2,3-DPG phosphatase deficiency  
Respiratory distress syndrome

notes

## disseminated intravascular coagulation screening (DIC screening)

**Type of test** Blood

### Normal findings

No evidence of DIC

### Test explanation and related physiology

This is a group of tests used to detect disseminated intravascular coagulation (DIC). Many pathologic conditions can instigate or are associated with DIC. The more common ones include bacterial septicemia, amniotic fluid embolism, retention of a dead fetus, malignant neoplasia, liver cirrhosis, extensive surgery (especially on the liver), post-extracorporeal heart bypass, extensive trauma, severe burns, and transfusion reactions.


In DIC, the entire clotting mechanism is inappropriately triggered. This results in significant systemic or localized intravascular formation of fibrin clots. The consequences of this futile clotting are small blood vessel occlusion and excessive bleeding caused by consumption of the platelets and clotting factors that have been used in intravascular clotting. The fibrinolytic system is also activated to break down the clot formation and the fibrin involved in the intravascular coagulation. This fibrinolysis results in the formation of fibrin degradation products (FDPs), which by themselves act as anticoagulants; these FDPs serve only to enhance the bleeding tendency.

Organ injury can occur as a result of the intravascular clots, which cause microvascular occlusion in various organs. This may cause serious anoxic injury in affected organs. Also, red blood cells passing through partly plugged vessels are injured and subsequently hemolyzed. The result may be ongoing hemolytic anemia. [Figure 15](#) summarizes DIC pathophysiology and effects.

When a patient with a bleeding tendency is suspected of having DIC, a series of readily performed laboratory tests should be conducted ([Table 15](#)). These tests are discussed separately in this book.

### Abnormal findings

Disseminated intravascular coagulation

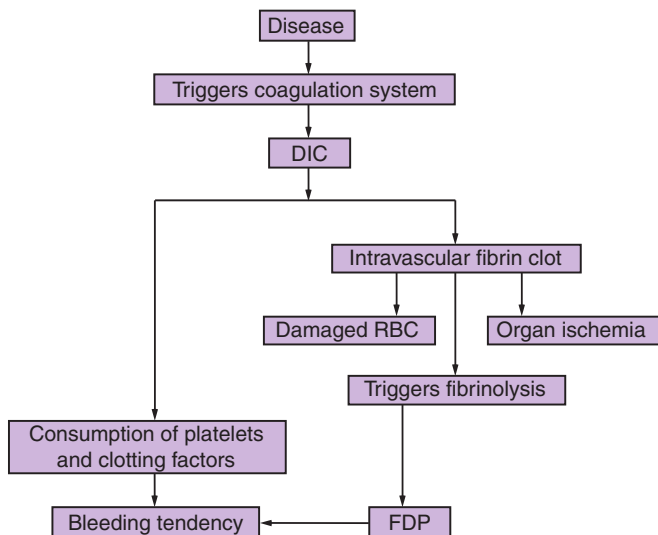

**FIGURE 15** Pathology of disseminated intravascular coagulation (DIC), which may result in bleeding tendency, organ ischemia, and hemolytic anemia. *FDP*, Fibrin degradation product; *RBC*, red blood cell.

**TABLE 15** Disseminated intravascular coagulation screening tests

Test	Result
Platelet count (see p. 706)	Decreased
Prothrombin time (see p. 753)	Prolonged
Partial thromboplastin time (see p. 681)	Prolonged
Coagulating factors (see p. 248)	Decreased factors I, II, V, VIII, X, and XIII (more commonly used for diagnosis rather than for screening)
Fibrin degradation products (see p. 877)	Increased
Fibrinogen (see p. 422)	Decreased
D-dimer (see p. 319)	Increased
Fibrinopeptide A (see p. 877)	Increased
Prothrombin fragment (see p. 877)	Increased

**drug monitoring** (Therapeutic drug monitoring [TDM], drug sensitivity genotyping)**Type of test** Blood**Normal findings**See [Table 16](#).**Test explanation and related physiology**

TDM entails measuring blood drug levels to determine effective drug dosages and to prevent toxicity. Drug monitoring is helpful in patients who take other medicines that may affect drug levels or act in a synergistic or antagonistic manner with the drug to be tested. There are some medicines (e.g., antiarrhythmics, bronchodilators, antibiotics, anticonvulsants, cardiotonics) that have a very narrow therapeutic margin (i.e., the difference between therapeutic and toxic drug levels is small).

[Table 16](#) (see p. 333) lists the therapeutic and toxic ranges for the average patient for some commonly tested drug levels. This list is far from complete. These ranges may not apply to all patients because clinical response is influenced by many factors. Also, note that different laboratories use different units for reporting test results and normal ranges. It is important that sufficient time pass between the administration of the medication and the collection of the blood sample to allow for adequate absorption and therapeutic levels to occur.

Blood is routinely used for TDM because results indicate what is presently going on with the drug at any one particular time. Urine drug levels reflect the presence of the drug over the previous several days.

Blood samples can be taken at the drug's *peak level* (highest concentration) or at the *trough level* (lowest concentration). Peak levels are useful when testing for toxicity, and trough levels are useful for demonstrating a satisfactory therapeutic level. Trough levels are often referred to as residual levels. The time when the sample should be drawn after the last dose of the medication varies according to whether a peak or trough level is requested and according to the half-life of the drug (the time required for the body to decrease the drug blood level by 50%). If peak levels are higher than the therapeutic range, toxicity may be experienced. If trough levels are lower than the therapeutic range, drug therapy is inadequate.



**TABLE 16 Drug monitoring data**

Drug	Use	Therapeutic level <sup>a</sup>	Toxic level <sup>a</sup>
Acetaminophen	Analgesic, antipyretic	Depends on use	> 25 mcg/mL
Aminophylline	Bronchodilator	10-20 mcg/mL	> 20 mcg/mL
Digoxin	Cardiac glycoside	0.8-2.0 ng/mL	> 2.4 ng/mL
Gentamicin	Antibiotic	5-10 mcg/mL	> 12 mcg/mL
Lidocaine	Antiarrhythmic	1.5-5.0 mcg/mL	> 5 mcg/mL
Lithium	Manic episodes of bipolar psychosis	0.8-1.2 mEq/L	> 2 mEq/L
Methotrexate	Antitumor agent	0.1-1 $\mu$ mol/24 hr	> 10 $\mu$ mol/24 hr
Phenobarbital	Anticonvulsant	10-30 mcg/mL	> 40 mcg/mL
Phenytoin	Anticonvulsant	10-20 mcg/mL	> 30 mcg/mL
Propranolol	Antiarrhythmic	50-100 ng/mL	> 150 ng/mL
Salicylate	Antipyretic, antiinflammatory, analgesic	100-250 mcg/mL	> 300 mcg/mL
Theophylline	Bronchodilator	10-20 mcg/mL	> 20 mcg/mL
Tobramycin	Antibiotic	5-10 mcg/mL	> 12 mcg/mL
Valproic acid	Anticonvulsant	50-100 mcg/mL	> 100 mcg/mL
Vancomycin	Antibiotic	Peak: 20-40 mcg/mL Trough: 5-15 mcg/mL	> 40 mcg/mL

<sup>a</sup> Levels vary according to the institution performing the test.

## Pharmacogenetics (genotyping for drug monitoring)

TDM is used to alter the dosage of medications to maximize efficacy and minimize side effects. There are several factors that affect efficacy and toxicity: patient compliance, patient age and size, access to adequate care, optimal dosing, and drug pharmacology issues (e.g., absorption, elimination, and drug interactions). Drugs undergo metabolism by enzyme systems to activate a bound (proactive) drug or to deactivate an active drug. The effectiveness of these enzyme systems of metabolism are determined by the genetic makeup of the patient. Based on the results of genotyping, a phenotype is predicted for a gene. Assessment of multiple genes (gene panels) may assist the clinician with personalized drug recommendations, avoidance of adverse drug reactions, and optimization of drug treatment.

Using the cytochromes as an example, four categories of drug metabolizers can be identified through genotype testing:

- Poor metabolizers (PMs)
- Intermediate metabolizers (IMs)
- Extensive metabolizers (EMs)
- Ultrarapid metabolizers (UMs)

Overall, whereas PMs and, to a lesser extent, IMs are prone to exaggerated side effects from active drugs, normal doses of the same drugs tend to be ineffectual for UMs. If a proactive drug is administered and must be hydrolyzed to its active form, PMs will not benefit from normal doses, but UMs will experience drug benefit from even small doses.

The cytochrome P450 (CYP450) system is a major family of drug-metabolizing enzymes. Several CYP450 enzymes are involved in the metabolism of a significant proportion of drugs. Drugs affected by cytochrome metabolism include medications used for allergy, pain, inflammatory, cardiovascular, cancer, viral, neurologic, and prostate diseases. Other medications using the cytochrome system for metabolism are anticoagulants, immunosuppressives, antacids, and sedative hypnotics.

Thiopurine methyltransferase (TPMT) and Nudix hydrolase (NUDT15) are other examples of metabolic enzyme systems that are specifically used in the metabolism of thiopurine drugs (e.g., azathioprine, 6-mercaptopurine [6MP], and 6-thioguanine).

Pharmacogenetics allows physicians to consider genetic information from patients in selecting medications and dosages of medications for a variety of common conditions (e.g., cardiac disease, psychiatric disease, and cancer).

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no food or fluid restrictions are needed.
- For patients suspected of having symptoms of drug toxicity, the best time to draw the blood specimen is when the symptoms are occurring.
  - If there is a concern regarding whether an adequate dose of the drug is achieved, it is best to obtain trough levels.

### During

- Collect a venous blood sample in a tube designated by the laboratory. Peak levels are usually obtained 1 to 2 hours after oral intake, approximately 1 hour after intramuscular (IM) administration, or approximately 30 minutes after intravenous (IV) administration. Residual (trough) levels are usually obtained shortly before (0-15 minutes) the next scheduled dose. Consult with the pharmacy for specific times.

### After

- Apply pressure or a pressure dressing to the venipuncture site.
- Clearly mark all blood samples with the following information: patient's name, diagnosis, name of drug, time of last drug ingestion, time of sample, and any other medications the patient is currently taking.
- Promptly send the specimen to the laboratory.

## Abnormal findings

Nontherapeutic levels of drugs

Toxic levels of drugs

notes

## ductoscopy (Mammary ductoscopy)

**Type of test** Endoscopy

### Normal findings

No tumor or premalignant changes

### Test explanation and related physiology

Most breast cancers start in the cells that line the milk ducts within the breast. Mammary ductoscopy refers to a procedure in which a miniaturized endoscope is used to get a closer look at the lining of the milk ducts and to provide access for biopsy or cell retrieval. Ductoscopy is used to visualize the breast ducts in women who have nipple discharge.

Breast diseases, including cancers, can be found at their very earliest stages with the use of this technique. Ductoscopy can identify cancers so small that mammography, ultrasound, and even MRI cannot see them. With this technique (Figure 16), premalignant changes can be identified and treated in an attempt to prevent breast cancer.

Breast ductal lavage (see p. 179) is a technique used to obtain and identify premalignant, atypical cells from breast ducts in patients who are considered at high risk for cancer and who have no evidence of breast malignancy on a mammogram or ultrasound. Ductoscopy is used in the hopes of identifying the causes of these changes (e.g., intraductal papillomas and early cancers) and possibly delivering ablative therapies to eradicate them.

### Interfering factors

- The inability to access the duct precludes performance of this endoscopic procedure.

### Procedure and patient care


#### Before

**PT** Explain the procedure to the patient.

- Be sure the breast examination and mammogram are normal.
- Obtain informed consent.

**PT** If the procedure is to be performed under general anesthesia, instruct the patient to abstain from eating and drinking for at least 8 hours.

- If the procedure is to be performed under local anesthesia, apply a topical anesthetic to the nipple area about ½ hour to 1 hour before the test.

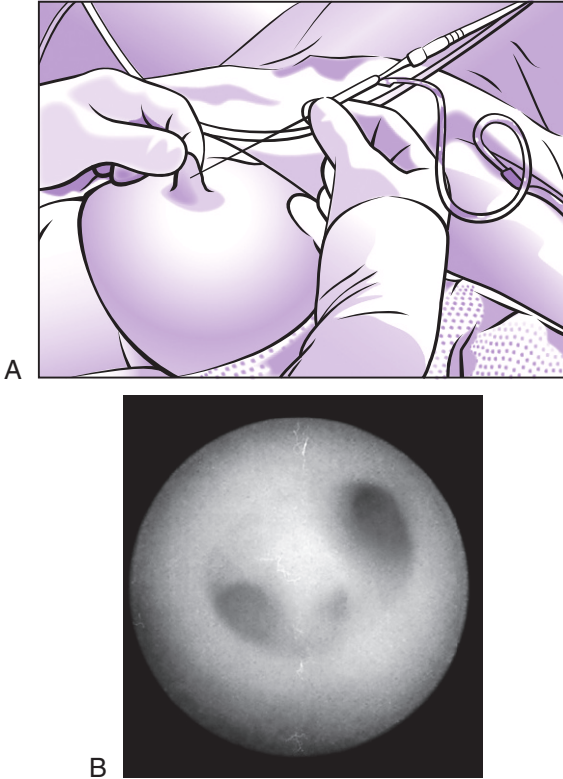

**FIGURE 16** A, Ductoscope is passed into the breast nipple. B, Image of normal ducts in the breast.

### During

- Note the following procedural steps:
  1. The breast is massaged to promote the discharge of nipple fluid. This helps to visually identify the ductal orifice in the nipple for endoscopy.
  2. The ductal opening in the nipple is gently dilated with progressively larger dilators. The mammary sheath containing the ductoscope is inserted and advanced under direct visualization as saline is injected to dilate the branches of the duct.
  3. The ductoscopy findings can be recorded on videotape.

### 338 ductoscopy

4. If any disease is identified, the scope can pinpoint the area for directed surgical removal.
  5. Ductal washings can also be obtained by aspirating some of the fluid for microscopic analysis.
- This procedure is usually performed by a surgeon in the office in approximately 30 minutes.

#### After

- PT** Inform the patient to contact the physician if she develops any redness, breast pain, or elevated temperature, which may indicate mastitis.

#### Abnormal findings

Invasive ductal cancer

Noninvasive ductal cancer

Atypical ductal hyperplasia

Papilloma

#### notes

## echocardiography (Cardiac echo, Heart sonogram, Transthoracic echocardiography [TTE])

**Type of test** Ultrasound

### Normal findings

Normal position, size, and movement of the cardiac valves and heart muscle wall

Normal directional flow of blood within the heart chambers

### Test explanation and related physiology

Echocardiography is a noninvasive ultrasound procedure used to evaluate the structure and function of the heart. *M-mode echocardiography* is a linear tracing of the motion of the heart structures over time. This allows the various cardiac structures to be located and studied regarding their movement during a cardiac cycle.

*Two-dimensional echocardiography* angles a beam within one sector of the heart. This produces an image of the spatial anatomic relationships within the heart.

*Three-dimensional echocardiography* allows for improved images of the heart wall and valves. The addition of high temporal resolution further improves images.

*Color Doppler echocardiography* detects the pattern of the blood flow and measures changes in velocity of blood flow within the heart and great vessels. Turbulent blood or altered velocity and direction of blood flow can be identified by changes in color. The most useful application of the color flow imaging is in determining the direction and turbulence of blood flow across regurgitant or narrowed valves. Doppler color flow imaging also may be helpful in assessing proper functioning of prosthetic valves.

Echocardiography, in general, is used in the diagnosis of a pericardial effusion, valvular heart disease (e.g., mitral valve prolapse, stenosis, regurgitation), subaortic stenosis, myocardial wall abnormalities (e.g., cardiomyopathy), infarction, and aneurysm. Cardiac tumors (e.g., myxomas) are easily diagnosed with ultrasound. Atrial and ventricular septal defects and other congenital heart diseases are also recognized by ultrasound. Finally, postinfarction mural thrombi are readily apparent with this testing.

Echocardiography is also used in *cardiac stress testing*. During an exercise or chemical cardiac stress test, ischemic muscle areas are evident as hypokinetic areas within the myocardium.

## 340 echocardiography

Echocardiography is being used with increasing frequency in emergent and urgent evaluations of patients with chest pains.

Echocardiography can be performed via the esophagus with use of a probe mounted on an endoscope. This is referred to as *transesophageal echocardiography* or *TEE* (see p. 903).

### Contraindications

- Patients who are uncooperative

### Interfering factors

- Chronic obstructive pulmonary disease (COPD)  
Patients who have severe COPD have a significant amount of air space between the heart and the chest cavity. Air space does not conduct ultrasound waves well.
- Obesity  
In obese patients, the space between the heart and the transducer is greatly enlarged; therefore accuracy of the test is decreased.

### Procedure and patient care

#### Before

**PT** Assure the patient that this is a painless study.

- Complete the request for the echocardiogram, including the pertinent patient history.

#### During

- Note the following procedural steps:
  1. The patient is placed in the supine position.
  2. Electrocardiographic (ECG) leads are placed (see p. 342).
  3. A gel, which allows better transmission of sound waves, is placed on the chest wall immediately under the transducer.
  4. Ultrasound is directed to the heart, and appropriate tracings are obtained.
- Note that this procedure usually takes approximately 45 minutes and is performed by an ultrasound technician in a darkened room within the cardiac laboratory or radiology department.

**PT** Tell the patient that no discomfort is associated with this study.

#### After

- Remove the gel from the patient's chest wall.

**PT** Inform the patient that the physician must interpret the study and that the results will be available in a few hours.



**Abnormal findings**

Valvular stenosis  
Valvular regurgitation  
Mitral valve prolapse  
Pericardial effusion  
Ventricular or atrial mural thrombi  
Myxoma  
Poor ventricular muscle motion  
Septal defects  
Ventricular hypertrophy  
Endocarditis

**notes**

**electrocardiography** (Electrocardiogram [ECG, EKG])

**Type of test** Electrodiagnostic

**Normal findings**


Normal rhythm, wave deflections, and heart rate (60-100 beats/min)

**Test explanation and related physiology**

The ECG is a graphic representation of the electrical impulses that the heart generates during the cardiac cycle. The monitoring electrodes detect the electrical activity of the heart from a variety of spatial perspectives. A *12-lead ECG* provides a comprehensive view of the flow of the heart's electrical currents in two different planes. There are six limb leads (combination of electrodes on the extremities) and six chest leads (corresponding to six sites on the chest). The limb leads provide a *frontal plane view* that bisects the body, separating the front and back; the chest leads provide a *horizontal plane view* that bisects the body, separating the top and bottom (Figure 17).


The normal ECG pattern is composed of waves arbitrarily designated by the letters *P*, *Q*, *R*, *S*, *T*, and *U*. The *Q*, *R*, and *S* waves are grouped together and described as the *QRS complex*. The significance of the waves and time intervals is as follows (Figure 18):

- *P wave*. This represents atrial electrical depolarization associated with atrial contraction. If *P* waves are absent or altered, the cardiac impulse originates outside the sinoatrial (SA) node.
- *PR interval*. This represents the time required for the impulse to travel from the SA node to the atrioventricular (AV) node. If this interval is prolonged, a conduction delay exists in the AV node (e.g., a first-degree heart block). If the *PR* interval is shortened, the impulse must have reached the ventricle through a “shortcut” (e.g., Wolff-Parkinson-White syndrome).
- *QRS complex*. This represents ventricular electrical depolarization associated with ventricular contraction. This complex consists of an initial downward (negative) deflection (*Q* wave), a large upward (positive) deflection (*R* wave), and a small downward deflection (*S* wave). A widened *QRS* complex indicates abnormal or prolonged ventricular depolarization time (e.g., a bundle-branch block), Wolff-Parkinson-White syndrome, or pacemaker rhythms.

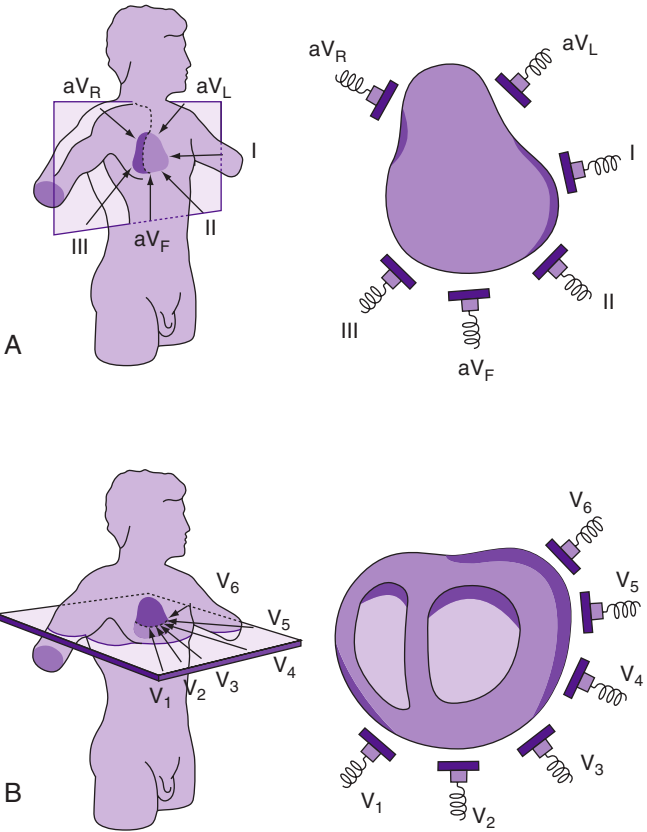

**FIG. 17** Planes of reference. **A**, The frontal plane. **B**, The horizontal plane.

- *ST segment*. This represents the period between the completion of depolarization and the beginning of repolarization of the ventricular muscle. This segment may be elevated or depressed in transient muscle ischemia (e.g., angina) or in muscle injury (e.g., the early stages of myocardial infarction).
- *T wave*. This represents ventricular repolarization (i.e., return to the resting state).

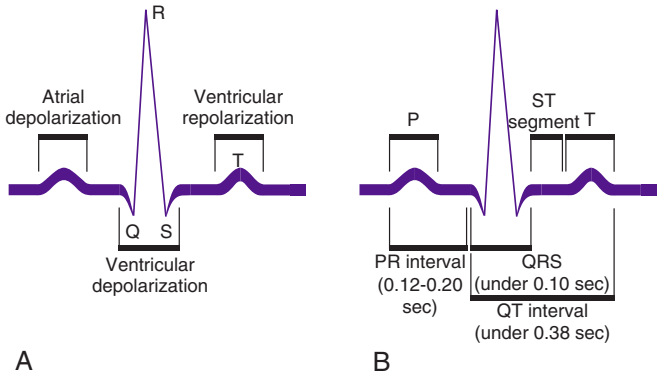

**FIG. 18** Electrocardiography. **A**, Normal electrocardiographic (ECG) deflections during depolarization and repolarization of the atria and ventricles. **B**, Principal ECG intervals between P, QRS, and T waves.

- *U wave.* This deflection follows the T wave and is usually quite small. It represents repolarization of the Purkinje nerve fibers within the ventricles.

The ECG is used primarily to identify abnormal heart rhythms (arrhythmias [dysrhythmias]) and to diagnose acute myocardial infarction, conduction defects, and ventricular hypertrophy. It is important to note that the ECG may be normal, even in the presence of heart disease.

For some patients at high risk for malignant ventricular arrhythmias, a *signal-averaged ECG (SAECG)* can be performed. This test averages several hundred QRS waveforms to detect late potentials that are likely to lead to ventricular arrhythmias. SAECGs have been a useful precursor to electrophysiologic (EPS) studies (see p. 359) because they can identify patients with unexplained syncope who may have ventricular tachycardias induced by the EPS study.

*Microvolt T-wave alternans (MTWA)* detects T-wave alternans (variations in the vector and amplitude of the T waves) on ECG signals as small as one-millionth of a volt. MTWA is defined as an alteration in the morphology of the T wave in an every-other-beat pattern. It has long been associated with ventricular arrhythmias and sudden death. MTWA is linked to the rapid onset of ventricular tachyarrhythmias.

MTWA is significant in the clinical context because it acts as a risk stratifier between patients who need implantable cardiac

defibrillators (ICDs) and those who do not. Patients who test negative for MTWA have a very low risk of sudden cardiac death and are less likely to require ICDs than those who test positive.

### Interfering factors

- Inaccurate placement of the electrodes
- Electrolyte imbalances
- Poor contact between the skin and the electrodes
- Movement or muscle twitching during the test
- ✚ Drugs that can affect results include barbiturates, digitalis, and quinidine.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no food or fluid restriction is necessary.
- PT** Assure the patient that the flow of electric current is *from* the patient. He or she will feel nothing during this procedure.
- Expose only the patient's chest and arms. Keep the abdomen and thighs adequately covered.

#### During

- Note the following procedural steps:
  1. The skin areas designated for electrode placement are prepared by using alcohol swabs or sandpaper to remove skin oil or debris. Sometimes the skin is shaved if the patient has a large amount of hair.
  2. Pads with special gel are applied to ensure electrical conduction between the skin and the electrodes.
  3. Electrodes are applied to the four extremities. Many cardiologists recommend that arm electrodes be placed on the upper arm because fewer muscle tremors are detected there.
  4. The chest leads are applied one at a time, three at a time, or six at a time, depending on the type of ECG machine. These leads are positioned as follows:
    - $V_1$ : in the fourth intercostal space (4ICS) at the right sternal border
    - $V_2$ : in 4ICS at the left sternal border
    - $V_3$ : midway between  $V_2$  and  $V_4$
    - $V_4$ : in 5ICS at the midclavicular line
    - $V_5$ : at the left anterior axillary line at the level of  $V_4$  horizontally
    - $V_6$ : at the left midaxillary line on the level of  $V_4$  horizontally

## 346 electrocardiography

- Note that cardiac technicians, nurses, or physicians perform this procedure in less than 5 minutes at the bedside or in the cardiology clinic.
- ▶ Tell the patient that although this procedure carries no discomfort, he or she must lie still in the supine position without talking while the ECG is recorded.

### After

- Remove the electrodes from the patient's skin and wipe off the electrode gel.
- Indicate on the ECG strip or request slip if the patient was experiencing chest pain during the study. The pain may be correlated with an arrhythmia on the ECG.

### Abnormal findings

Cardiac arrhythmias

Acute myocardial infarction

Myocardial ischemia

Old myocardial infarction

Conduction defects

Conduction system disease

Wolff-Parkinson-White syndrome

Ventricular hypertrophy

Cor pulmonale

Pulmonary embolus

Electrolyte imbalance

Pericarditis

### notes

## electroencephalography (Electroencephalogram [EEG])

**Type of test** Electrodiagnostic

### Normal findings

Normal frequency, amplitude, and characteristics of brain waves

### Test explanation and related physiology

The EEG is a graphic recording of the electrical activity of the brain. EEG electrodes are placed on the scalp over multiple areas of the brain to detect and record electrical impulses within the brain. This study is invaluable in the investigation of epileptic states to detect the focus of seizure activity. Patients with cerebral lesions (e.g., tumors, infarctions) have abnormally slow EEG waves, depending on the size and location of the lesion. Because this study determines the overall activity of the brain, it can be used to evaluate trauma and drug intoxication and to determine cerebral death in comatose patients.

The EEG also can be used to monitor the EPS effects of cerebral blood flow. For example, during carotid endarterectomy, the carotid vessel must be temporarily occluded. When this surgery is performed with the patient under general anesthesia, the EEG can be used for early detection of cerebral tissue ischemia. Temporary shunting of the blood during the surgery is then required.

*Electrocorticography (ECoG)* is a form of EEG performed during craniotomy. Electrodes are placed directly on the exposed surface of the brain to record electrical activity from the cerebral cortex. ECoG is currently considered to be the standard for defining epileptogenic zones before attempts at surgical interruption are made. This procedure is invasive. The same information can be obtained by a noninvasive brain imaging technique called *magnetoencephalography (MEG)*.

MEG is a noninvasive imaging technique used to measure the magnetic fields produced by electrical activity in the brain with an extremely sensitive device called a *superconducting quantum interference device (SQUID)*. The data obtained by MEG are commonly used to assist neurosurgeons in localizing pathology or defining sites of origin for epileptic seizures. MEG is also used in localizing important adjacent cortical areas for surgical planning in patients with brain tumors or intractable epilepsy.

### Interfering factors

- Fasting may cause hypoglycemia, which could modify the EEG pattern.
- Drinks containing caffeine (e.g., coffee, tea, cocoa, cola) interfere with test results.
- Body and eye movements during the test can cause changes in brain wave patterns.
- ✔ Drugs that may affect test results include sedatives.

### Procedure and patient care

#### Before

- PT Explain the procedure to the patient.
- PT Assure the patient that this test cannot read the mind or detect senility.
- PT Assure the patient that the flow of electrical activity is *from* the patient. He or she will not feel anything.
- PT Instruct the patient to wash his or her hair the night before the test. No oils, sprays, or lotion should be used.
  - Check to see if the physician wants to discontinue any medications before the study.
- PT Instruct the patient if sleeping time should be shortened the night before the test. Adults may not be allowed to sleep more than 4 or 5 hours and children not more than 5 to 7 hours if a sleep EEG will be done.

#### During

- Note the following procedural steps:
  1. The EEG is usually performed in a specially constructed room that is shielded from outside disturbances.
  2. The patient is placed in a supine position on a bed or reclining in a chair.
  3. Sixteen or more electrodes are applied to the scalp with electrode paste in a uniform pattern over both sides of the head, covering the prefrontal, frontal, temporal, parietal, and occipital areas.
  4. One electrode may be applied to each earlobe for grounding.
  5. After the electrodes are applied, the patient is instructed to lie still with eyes closed.
  6. The technician observes the patient during the EEG recording for movements that could alter results.
  7. Approximately every 5 minutes, the recording is interrupted to permit the patient to move if desired.



- In addition to the resting EEG, note that the following *activating procedures* can be performed:
    1. The patient is *hyperventilated* (asked to breathe deeply 20 times a minute for 3 minutes) to induce alkalosis and cerebral vasoconstriction, which can activate abnormalities.
    2. *Photostimulation* is performed by flashing a light over the patient's face with the eyes opened or closed. Seizure activity may be seen on the EEG.
    3. A *sleep EEG* may be performed to aid in the detection of some abnormal brain waves that are seen only if the patient is sleeping (e.g., frontal lobe epilepsy). The sleep EEG is performed after orally administering a sedative-hypnotic. A recording is performed when the patient is falling asleep, when the patient is sleeping, and when the patient is awakening.
  - Note that this study is performed by an EEG technician in approximately 45 minutes to 2 hours.
- PT** Tell the patient that no discomfort is associated with this study, other than possibly missing sleep.

### After

- Help the patient remove the electrode paste.
- Ensure safety precautions until the effects of any sedatives have worn off. Keep the bed's side rails up.

### Abnormal findings

Seizure disorders (e.g., epilepsy)

Brain tumor

Brain abscess

Head injury

Cerebral death

Encephalitis

Intracranial hemorrhage

Cerebral infarct

Narcolepsy

Alzheimer disease

### notes

**electromyography (EMG)**

**Type of test** Electrodiagnostic

**Normal findings**

No evidence of neuromuscular abnormalities

**Test explanation and related physiology**

EMG is used to detect primary muscular disorders, along with muscular abnormalities caused by other system diseases (e.g., nerve dysfunction, sarcoidosis, paraneoplastic syndrome).

Spontaneous muscle movement such as fibrillation and fasciculation can be detected during EMG. When evident, these waveforms indicate injury or disease of the nerve or muscle being evaluated. A decrease in the number of muscle fibers able to contract is typically observed with peripheral nerve damage. This study is usually done in conjunction with nerve conduction studies (see p. 355) and may also be called electromyoneurography.

**Contraindications**

- Patients receiving anticoagulant therapy
- Patients with extensive skin infection

**Potential complications**

- Rarely, hematoma at the needle insertion site

**Interfering factors**

- Edema, hemorrhage, or thick subcutaneous fat can interfere with test results.
- Patients with excessive pain may have false results.

**Procedure and patient care****Before**

**PT** Explain the procedure to the patient. Allay any fears and allow the patient to express concerns.

- Obtain informed consent if required by the institution.

**PT** Tell the patient that fasting is not usually required; however, some facilities may restrict stimulants (coffee, tea, cocoa, cola, cigarettes) for 2 to 3 hours before the test.

- If serum enzyme tests (e.g., aspartate transaminase, creatinine phosphokinase, lactate dehydrogenase) are ordered, the specimen should be drawn before EMG or 5 to 10 days after the test because the EMG may cause misleading elevations of these enzymes.

- Premedication or sedation is usually avoided because of the need for patient cooperation.

### During

- Note the following procedural steps:
    1. This study is usually done in an EMG laboratory.
    2. The patient's position depends on the muscle being studied.
    3. A tiny needle that acts as a recording electrode is inserted into the muscle being examined or overlying the nerve itself. In most circumstances, however, that reference electrode is in the needle itself.
    4. A reference electrode is placed nearby on the skin surface.
    5. The patient is asked to keep the muscle at rest.
    6. The oscilloscope display is viewed for any evidence of spontaneous electrical activity, such as fasciculation or fibrillation.
    7. The patient is asked to contract the muscle slowly and progressively.
    8. The electrical waves produced are examined for their number, form, and amplitude.
    9. A nerve innervating a particular muscle group is stimulated, and the resulting muscle contraction is evaluated as described if nerve conduction studies are performed concomitantly.
  - Note that the EMG is performed by a physical therapist, physiatrist, or neurologist in approximately 20 minutes.
- PT** Tell the patient that the small needle size makes this procedure nearly painless.

### After

- Observe the needle site for hematoma or inflammation.
- Provide pain medication if needed.

**Abnormal findings**

Polymyositis

Muscular dystrophy

Myopathy

Traumatic injury

Hyperadrenalism

Hypothyroidism

Paraneoplastic syndrome (e.g., lung cancer)

Sarcoidosis

Guillain-Barré syndrome

Myasthenia gravis

Peripheral nerve injury, entrapment, or compression

Spinal cord injury or disease

Acetylcholine blockers (e.g., curare, snake venom)

Multiple sclerosis

Diabetic neuropathy

Anterior poliomyelitis

Muscle denervation

Amyotrophic lateral sclerosis

**notes**

## electromyography of the pelvic floor sphincter

(Pelvic floor sphincter electromyography, Pelvic floor sphincter EMG, Rectal EMG procedure)

**Type of test** Electrodiagnostic

### Normal findings

Increased EMG signal during bladder filling

Silent EMG signal on voluntary micturition

Increased EMG signal at the end of voiding

Increased EMG signal with voluntary contraction of the anal sphincter

### Test explanation and related physiology

This test uses the placement of electrodes on or in the pelvic floor musculature to evaluate the neuromuscular function of the urinary or anal sphincter. It is performed most often in patients who have urinary or fecal incontinence. The pathology causing the muscle weakness can be muscular or neurologic. With pelvic floor sphincter EMG, these two causes can be separated.

The main benefit of this study is to evaluate the external sphincter (skeletal muscle) activity during voiding. This test is also used to evaluate the bulbocavernosus reflex and voluntary control of external sphincter or pelvic floor muscles. The pelvic floor sphincter EMG also aids in the investigation of functional or psychological disturbances of voiding. Fecal incontinence caused by muscular dysfunction can also be evaluated by rectal sphincter EMG.

Recordings may be made from surface or needle electrodes within the muscle; surface electrodes are most often used. These electrodes allow for observation of and change in the muscle activity before and during voiding.

### Contraindications

- Patients who cannot cooperate during the procedure

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

**PT** Inform the patient that cooperation is essential.

## 354 electromyography of the pelvic floor sphincter

### During

- Note the following procedural steps:
  1. Electrodes are placed at the 2 o'clock and 10 o'clock positions on the perianal skin to monitor the pelvic floor musculature during voiding.
  2. A third electrode is usually placed on the thigh and serves as a ground.
  3. Electrical activity is recorded with the bladder empty and the patient relaxed.
  4. Reflex activity is evaluated by asking the patient to cough and by stimulating the urethra and trigone by gently tugging on an inserted Foley catheter (bulbocavernous reflex).
  5. Voluntary activity is evaluated by asking the patient to contract and relax the sphincter muscle.
  6. The bladder is filled with sterile water at room temperature at a rate of 100 mL/min.
  7. The EMG responses to filling and detrusor hyperreflexia (if present) are recorded.
  8. Finally, when the bladder is full and with the patient in a voiding position, the filling catheter is removed and the patient is asked to urinate. In the normal patient, the EMG signals build during bladder filling and cease promptly on voluntary micturition, remaining silent until the pelvic floor contracts at the end of voiding.
  9. The electrical waves produced are examined for their number and form.
- Note that a urologist performs this study in less than 30 minutes.
- PT** Explain to the patient that this study is slightly more uncomfortable than urethral catheterization.

### After

- If needle electrodes were used, observe the needle sites for hematoma or inflammation.

### Abnormal findings

Neuromuscular dysfunction of lower urinary sphincter

Pelvic floor muscle dysfunction of anal sphincter

notes

**electroneurography** (ENG, Nerve conduction studies)**Type of test** Electrodiagnostic**Normal findings**

No evidence of peripheral nerve injury or disease

**Test explanation and related physiology**

Nerve conduction studies (NCS) evaluate the integrity of the nerves and allow for the detection and location of peripheral nerve injury or disease. By initiating an electrical impulse at one site (proximal when evaluating motor nerves or distal when evaluating sensory nerves) of a nerve and recording the time required for that impulse to travel to a second site (opposite above) of the same nerve, the conduction velocity of an impulse in that nerve can be determined. This study is usually done in conjunction with EMG (see p. 350) and also may be called *electromyoneurography*.

Because conduction velocity may require contraction of a muscle as an indication of an impulse arriving at the recording electrode, significant primary muscular disorders may cause a falsely slow nerve conduction velocity. This “muscular” variable is eliminated if one evaluates the suspected pathologic muscle group before performing NCS. This muscular factor is evaluated by measuring distal latency (i.e., the time required for stimulation of the nerve to cause muscular contraction). As the distal latency is calculated, the motor NCS is performed normally by stimulating the nerve bundle. Conduction velocity can then be determined by the following equation:

$$\text{Conduction velocity (in meters per second)} = \frac{\text{Distance (in meters)}}{\text{Total latency} - \text{Distal latency}}$$

**Interfering factors**

- Patients in severe pain may have false results.

**Procedure and patient care****Before****PT** Explain the procedure to the patient. Allay any fears.

- Obtain informed consent if required by the institution.

**PT** Tell the patient that no fasting or sedation is usually required.**During**

- Note the following procedural steps:
  1. This test can be performed in a nerve conduction laboratory or at the patient’s bedside.

## 356 electroneurography

2. The patient's position depends on the area of suspected peripheral nerve injury or disease.
  3. A recording electrode is placed on the skin overlying a muscle innervated solely by the relevant nerve or overlying the nerve itself. All skin-to-electrode connections are ensured by using electrical gel.
  4. A reference electrode is placed nearby.
  5. The nerve is stimulated by a shock-emitting device at an adjacent location.
  6. The time between nerve impulse and muscular contraction (distal latency) is measured in milliseconds.
  7. The nerve is similarly stimulated at a location proximal to the area of suspected injury or disease.
  8. The time required for the impulse to travel from the site of initiation to muscle contraction (total latency) is recorded in milliseconds.
  9. The distance between the site of stimulation and the recording electrode is measured in centimeters.
  10. Conduction velocity is converted to meters per second and is computed as in the previous equation.
- Note that this test takes approximately 40 minutes and is performed by a physiatrist or a neurologist.
- PT** Tell the patient that this test may be uncomfortable because a mild shock is required for nerve impulse stimulation.

### After

- Remove the electrode gel from the patient's skin.

### Abnormal findings

Peripheral nerve injury or disease

Carpal tunnel syndrome

Herniated disc disease

Poliomyelitis

Diabetic neuropathy

Myasthenia gravis

Guillain-Barré syndrome

### notes



**electronystagmography** (Electrooculography)

**Type of test** Electrodiagnostic

**Normal findings**

Normal nystagmus response

Normal oculovestibular reflex

**Test explanation and related physiology**

Electronystagmography is used to evaluate nystagmus (involuntary rapid eye movement) and the muscles controlling eye movement. By measuring changes in the electrical field around the eye, this study can make a permanent recording of eye movement at rest, with a change in head position, and in response to various stimuli. The test delineates the presence or absence of nystagmus, which is caused by the initiation of the oculovestibular reflex. Nystagmus should occur when initiated by positional, visual, or caloric stimuli. Unlike caloric studies (see p. 192), in which nystagmus is usually determined visually, electronystagmography electrically records the direction, velocity, and degree of nystagmus.

If nystagmus does not occur with stimulation, the vestibular cochlear apparatus, cerebral cortex (temporal lobe), auditory nerve, or brainstem is abnormal. Tumors, infection, ischemia, and degeneration can cause such abnormalities. When put together with the entire clinical picture, the pattern of nystagmus helps in the differentiation between central and peripheral vertigo.

This study also may help evaluate unilateral hearing loss and vertigo. Unilateral hearing loss may be caused by middle ear problems or nerve injury. If the patient experiences nystagmus with stimulation, the auditory nerve is working and hearing loss can be blamed on the middle ear.

**Contraindications**

- Patients with perforated eardrums, who should not have water irrigation
- Patients with pacemakers

**Interfering factors**

- Blinking of the eyes can alter test results.
- ✚ Drugs that can alter results include antivertigo agents, sedatives, and stimulants.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient not to apply facial makeup before the test because electrodes will be taped to the skin around the eyes.
- PT** Instruct the patient not to eat solid food before the test to reduce the likelihood of vomiting.
- PT** Instruct the patient not to drink caffeine or alcoholic beverages for approximately 24 to 48 hours (as ordered).
- Check with the physician regarding withholding any medications that could interfere with the test results.

### During

- Note the following procedural steps:
  1. This procedure is usually performed in a darkened room with the patient seated or lying down.
  2. If there is any wax in the ear, it is removed.
  3. Electrodes are taped to the skin around the eyes.
  4. Various procedures are used to stimulate nystagmus such as pendulum tracking, changing head position, changing gaze position, and caloric tests (see p. 192).
  5. Several recordings are made with the patient at rest and to demonstrate patient response to various procedures (e.g., blowing air into the ear, irrigating the ear with water).
  6. Nystagmus response is compared with the expected ranges, and the results are recorded as normal, borderline, or abnormal.
- Note that this procedure is performed by a physician or audiologist in approximately 1 hour.
- PT** Tell the patient that nausea and vomiting may occur.

### After

- Consider prescribing bed rest until nausea, vertigo, or weakness subsides.

## Abnormal findings

Brainstem lesions  
Cerebellum lesions  
Auditory nerve damage  
Vestibular system lesions  
Congenital disorder  
Demyelinating disease

### notes

**electrophysiologic study** (EPS, Cardiac mapping)

**Type of test** Electrodiagnostic; manometric

**Normal findings**

Normal conduction intervals, refractive periods, and recovery times

**Tilt-table testing**

< 20 mm Hg decrease in systolic blood pressure

< 10 mm Hg increase in diastolic blood pressure

Heart rate increase should be less than 10 beats/min.

**Test explanation and related physiology**

In this invasive procedure, multiple electrode catheters are fluoroscopically placed through a peripheral vein and into the right atrium and/or ventricle or, less often, through an artery into the left atrium and/or ventricle. With close cardiac monitoring, the electrode catheters are used to pace the heart and potentially induce arrhythmias. Defects in the heart conduction system can then be identified; arrhythmias that are otherwise not apparent also can be induced, identified, and treated. The effectiveness of antiarrhythmic drugs can be assessed.

EPS can also be therapeutic. With the use of radiofrequency, sites with documented low thresholds for inducing arrhythmias can be obliterated to stop the arrhythmias.

The *tilt-table test* is sometimes performed with an EPS cardiac study and is a provocative test used to diagnose vasopressor syncope syndrome. Patients with this syndrome usually demonstrate symptomatic hypotension and syncope within a few to 30 minutes of being tilted upright by approximately 60 to 80 degrees. The tilt-table test is often used to assess the efficacy of prophylactic pacing in some patients with vasopressor syncope. It is also used to evaluate the influence of posture with some forms of tachyarrhythmias. Normally, a minimal drop in systolic blood pressure, rise in diastolic blood pressure, and increase in heart rate occur in the tilted position. Patients with vasopressor syncope demonstrate these changes in an exaggerated fashion and become lightheaded and dizzy on assuming the tilted position.

**Contraindications**

- Patients who are uncooperative
- Patients with acute myocardial infarction

## 360 electrophysiologic study

### Potential complications

- Cardiac arrhythmias leading to ventricular tachycardia or fibrillation
- Perforation of the myocardium
- Catheter-induced embolic cerebrovascular accident (CVA; stroke) or myocardial infarction
- Peripheral vascular problems
- Hemorrhage
- Phlebitis at the venipuncture site

### Interfering factors

- Patients with dehydration or hypovolemia demonstrate similar changes in blood pressure and heart rate with tilt-table testing. This is especially true in elderly patients.
- ✚ Patients on antihypertensive medications or diuretics also may demonstrate similar changes when placed in the tilt position.
- ✚ Drugs that may interfere with test results include analgesics, sedatives, and tranquilizers.

### Procedure and patient care

#### Before

- PT** Instruct the patient to fast for 6 to 8 hours before the procedure. Usually fluids are permitted until 3 hours before.
- Obtain an informed consent from the patient.
- PT** Encourage the patient to verbalize any fears.
- Prepare the catheter insertion site as directed.
- Collect a blood sample for potassium or drug levels if indicated.
- Obtain peripheral IV access for the administration of drugs.
- Inquire as to whether the patient has had excessive fluid loss (diarrhea or vomiting) in the previous 24 hours.
- Record the use of antihypertensive or diuretic medicines.

#### During

- Note the following procedural steps:

#### *Electrophysiologic study*

1. In the cardiac catheterization laboratory, the patient has ECG leads attached.
2. The catheter insertion site, usually the femoral vein, is prepared and draped in a sterile manner.
3. Under fluoroscopic guidance, the catheter is passed to the atrium and ventricle.
4. Baseline surface intracardiac ECGs are recorded.
5. Various parts of the cardiac electroconduction system are stimulated by atrial or ventricular pacing.

6. Mapping of the electroconduction system is performed.
7. Arrhythmias are identified.
8. Drugs may be administered to assess their efficacy in preventing EPS-induced arrhythmias.
9. Because dangerous arrhythmias can be prolonged, cardioversion must be immediately available.
10. The patient is also constantly engaged in light conversation to assess mental status and consciousness.

#### *Tilt-table test*

1. The patient lies supine on a horizontal tilt table.
2. Obtain the patient's blood pressure and pulse as baseline values before tilting is carried out.
3. Monitor these vital signs during the procedure.
4. Question the patient as to the presence of symptoms of dizziness and lightheadedness.
5. The table is progressively tilted 60 to 80 degrees while the patient is being monitored. Alternatively, the patient is asked to sit or stand.

**PT** Tell the patient that he or she may experience palpitations, lightheadedness, or dizziness when arrhythmias are induced. For most patients, this is an anxiety-producing experience.

**PT** Inform the patient that discomfort from catheter insertion is minimal.

#### **After**

- Keep the patient on bed rest for approximately 6 to 8 hours.
- Evaluate the access site for swelling and bleeding.
- Monitor the patient's vital signs for at least 2 to 4 hours for hypotension and arrhythmias.
- Continue cardiac monitoring to identify arrhythmias.
- Cover the area with sterile dressings if the electrical catheter is left in place for subsequent studies.

#### **Abnormal findings**

Electroconduction defects

Cardiac arrhythmias

SA node defects (e.g., sick sinus syndrome)

Atrioventricular node defects and heart blocks

Vasomotor syncope syndrome

Inducible arrhythmias (e.g., ventricular tachycardia and Wolff-Parkinson-White syndrome)

## endometrial biopsy

**Type of test** Microscopic examination of tissue

### Normal findings

No pathologic conditions

Presence of a secretory-type endometrium 3 to 5 days before normal menses

### Test explanation and related physiology

An endometrial biopsy can determine whether ovulation has occurred. A biopsy specimen taken 3 to 5 days before normal menses should demonstrate a *secretory-type* endometrium on histologic examination if ovulation and corpus luteum formation have occurred. If not, only a preovulatory *proliferative-type* endometrium will be seen.

This test can determine whether a woman has adequate ovarian estrogen and progesterone levels. Another major use of endometrial biopsy is to diagnose endometrial cancer, tuberculosis, polyps, or inflammatory conditions and to evaluate dysfunctional uterine bleeding.

### Contraindications

- Patients with infections (e.g., trichomonal, candidal, or suspected gonococcal) of the cervix or vagina
- Patients in whom the cervix cannot be visualized (e.g., because of abnormal position or previous surgery)
- Patients who are or may be pregnant because the procedure may induce labor or abortion

### Potential complications

- Perforation of the uterus
- Uterine bleeding
- Interference with early pregnancy
- Infection

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Ensure that written and informed consent is obtained.

**PT** Tell the patient that neither fasting nor sedation is usually required.

## During

- Note the following procedural steps:
    1. The patient is placed in the lithotomy position, and a pelvic examination is performed to determine the position of the uterus.
    2. The cervix is exposed and cleansed.
    3. A biopsy instrument is inserted into the uterus, and specimens are obtained from the anterior, posterior, and lateral walls. This can be performed as part of a dilation and curettage or hysteroscopy (see p. 533).
    4. The specimens are placed in a solution containing 10% formalin solution and sent to pathology.
  - Note that this procedure is performed by an obstetrician/gynecologist in approximately 10 to 30 minutes.
- PT** Tell the patient that this procedure may cause momentary discomfort (menstrual-type cramping).

## After

- Assess the patient's vital signs at regular intervals. Any temperature elevation should be reported to the physician because this procedure may activate pelvic inflammatory disease.
- PT** Advise the patient to wear a pad because some vaginal bleeding is to be expected. Tell the patient to call her physician if there is excessive bleeding (> 1 pad per hour).
- PT** Inform the patient that douching and intercourse are not permitted for 72 hours after the biopsy.
- PT** Instruct the patient to rest during the next 24 hours and to avoid heavy lifting to prevent uterine hemorrhage.

## Abnormal findings

Anovulation

Tumor

Tuberculosis

Polyps

Inflammatory condition

## notes

## endoscopic retrograde cholangiopancreatography (ERCP)

**Type of test** Endoscopy

### Normal findings

Normal size of biliary and pancreatic ducts

No obstruction or filling defects in the biliary or pancreatic ducts

### Test explanation and related physiology

With the use of a fiberoptic endoscope, ERCP provides radiographic visualization of the biliary and pancreatic ducts. This is especially useful in patients with jaundice. If a partial or total obstruction of these ducts exists, characteristics of the obstructing lesion can be demonstrated. Stones, benign strictures, cysts, ampullary stenosis, anatomic variations, and malignant tumors can be identified.

Incision of the papillary muscle in the ampulla of Vater can be performed through the scope at the time of ERCP. This incision widens the distal common duct so that common bile duct gallstones can be removed. Stents can be placed through narrowed bile ducts with the use of ERCP, and the bile of jaundiced patients can be internally drained. Pieces of tissue and brushings of the common bile duct can be obtained by ERCP for pathologic review. Manometric studies of the sphincter of Oddi and pancreaticobiliary ducts can be performed at the time of ERCP. These are used to investigate unusual functional abnormalities of these structures.

Another less commonly used method of visualizing the biliary tree is *percutaneous transhepatic cholangiography (PTHC)*. PTHC is performed by passing a needle through the skin into the liver and into an intrahepatic bile duct. Iodinated x-ray contrast dye is directly injected into the biliary system. The intrahepatic and extrahepatic biliary ducts and occasionally the gallbladder can be visualized. If the jaundice is found to result from extrahepatic obstruction, a catheter can be left in the bile duct and used for external drainage of bile. Furthermore, with the assistance of ERCP, a stent can be placed across a stricture to decompress the biliary system internally.

Magnetic resonance imaging cholangiopancreatography (see p. 601) is a noninvasive test that can provide information similar to ERCP. However, interventional procedures such as papillotomy cannot be performed with this testing method.



## Contraindications

- Patients who are uncooperative. Cannulation of the ampulla of Vater requires that the patient lie very still.
- Patients whose ampulla of Vater is not accessible endoscopically because of previous upper gastrointestinal (GI) surgery
- Patients with esophageal diverticula. The scope can fall into a diverticulum and perforate its wall.
- Patients with known acute pancreatitis

## Potential complications

- Perforation of the esophagus, stomach, or duodenum
- Gram-negative sepsis. This results from introducing bacteria through the biliary system and into the blood.
- Pancreatitis resulting from pressure of the dye injection
- Aspiration of gastric contents into the lungs
- Respiratory arrest as a result of oversedation

## Interfering factors

- Barium in the abdomen as a result of a previous upper GI series or barium enema x-ray study precludes adequate visualization of the biliary and pancreatic ducts.

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.


- Obtain informed consent from the patient.

**PT** Inform the patient that breathing will not be compromised by the insertion of the endoscope.

- Keep the patient NPO as of midnight the day of the test.
- Administer appropriate premedication if ordered.

### During

- Note the following procedural steps:
  1. A flat plate of the abdomen is taken to ensure that any barium from previous studies will not obscure visualization of the bile duct.
  2. The patient is positioned supine or on the left side.
  3. The patient is usually sedated with a narcotic and a sedative-hypnotic.
  4. The pharynx is sprayed with a local anesthetic to inactivate the gag reflex and to lessen the discomfort.
  5. A fiberoptic duodenoscope is inserted through the oral pharynx and passed through the esophagus and stomach and then into the duodenum (Figure 19).

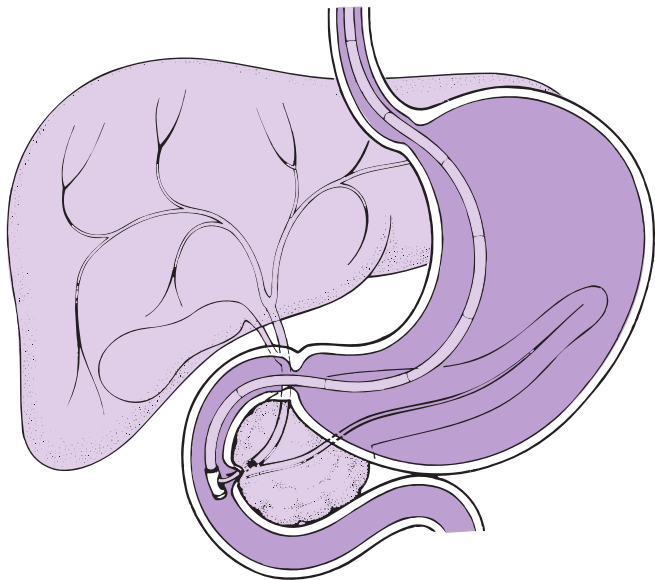

**FIG. 19** Endoscopic retrograde cholangiopancreatography. The fiberoptic scope is passed into the duodenum. Note the small catheter being advanced into the biliary duct.

6. Glucagon is often administered intravenously to minimize the spasm of the duodenum and to improve visualization of the ampulla of Vater.
  7. Through the accessory lumen within the scope, a small catheter is passed through the ampulla of Vater and into the common bile or pancreatic ducts.
  8. Radiographic dye is injected, and x-ray images are taken.
- Note that the test usually takes approximately 1 hour and is performed by a physician trained in endoscopy. The x-ray images are interpreted by a radiologist.
- PT** Tell the patient that no discomfort is associated with the dye injection but that minimal gagging may occur during the initial introduction of the scope into the oral pharynx.

#### After

- Do *not* allow the patient to eat or drink until the gag reflex returns. Encourage light eating for the next 12 to 24 hours.

- Observe the patient closely for development of abdominal pain, nausea, and vomiting. This may herald the onset of ERCP-induced pancreatitis or gastroduodenal perforation.
  - Observe safety precautions until the effects of the sedatives have worn off.
  - Monitor the patient for signs of respiratory depression. Medication (e.g., naloxone) should be available to counteract serious respiratory depression.
  - Assess the patient for signs and symptoms of septicemia, which may indicate the onset of ERCP-induced cholangitis.
- PT** Inform the patient that he or she may be hoarse and have a sore throat for several days. Drinking cool fluids and gargling will help relieve some of this soreness.
- PT** Instruct the patient to notify the physician immediately of fever or shaking chills. This may indicate possible cholangitis.

### Abnormal findings

Tumor, strictures, or gallstones of the common bile duct

Sclerosing cholangitis

Biliary sclerosis

Cysts of the common bile duct

Tumor, strictures, or inflammation of the pancreatic duct

Pseudocyst of the pancreatic duct

Chronic pancreatitis

Anatomic biliary or pancreatic duct variations

Cancer of the duodenum or ampulla of Vater

### notes

**Epstein-Barr virus testing (EBV antibody titer)**

**Type of test** Blood

**Normal findings**

Titers  $\leq 1:10$  are nondiagnostic.

Titers of 1:10 to 1:60 indicate infection at some undetermined time.

Titers of  $\geq 1:320$  suggest active infection.

Four-fold increase in titer in paired sera drawn 10 to 14 days apart is usually indicative of an acute infection.

**Test explanation and related physiology**

EBV infects 80% of the U.S. population. After infection occurs, the virus becomes latent but can be reactivated later. EBV infection can produce infectious mononucleosis. Mononucleosis is seen most often in children, adolescents, and young adults. Laboratory findings of lymphocytosis, atypical lymphocytes, and transient serum heterophil antibodies are seen in patients with acute EBV infection. In Africa, EBV has been associated with Burkitt lymphoma. In China, EBV infection has been associated with nasopharyngeal carcinoma.

The majority of EBV infections can be recognized by testing the patient's serum for heterophile antibodies (rapid latex slide agglutination test or mononucleosis [mono] spot test; see p. 630). Other more specific immunologic tests are recommended only when a mononucleosis screening procedure is negative and infectious mononucleosis or a complication of EBV infection is suspected. Also, more specific tests can more precisely define the acuity of the infection (Table 17). When EBV is suspected but the heterophile antibody is not detected, an evaluation of the EBV-specific antibody profile (e.g., EBV viral capsid antigen [VCA] IgM, EBV VCA IgG, and EBV nuclear antigen [EBNA]) may be useful (Table 18). The viral capsid antigen-antibodies (VCAs) can be immunoglobulin (Ig) G or IgM. The EBV nuclear antigen (EBNA) is located in the nuclei of the infected lymphocyte. Another EBV antigen is called the early antigen (EA). There are two EA antigens. One is EA-D and is commonly associated with nasopharyngeal cancer. EA-R is commonly associated with Burkitt lymphoma.

The interpretation of EBV antibody tests is based on these assumptions:

1. After the person becomes infected with EBV, the anti-VCA IgG antibodies appear first.

**TABLE 17** Serologic studies and the timing of infections

Serologic study	Appears/disappears	Clinical significance
Monospot heterophil	5 days/2 weeks	Acute or convalescent infection
VCA-IgM	7 days/3 months	Acute or convalescent infection
VCA-IgG	7 days/exists for life	Acute, convalescent, or old infection
EBNA-IgG	3 weeks/exists for life	Old infection
EA-D	7 days/2 weeks	Acute or convalescent infection

*EA*, Early antigen; *EA-D*, early antigen D; *EBNA*, Epstein-Barr virus nuclear antigen; *IgG*, immunoglobulin G; *IgM*, immunoglobulin M; *VCA*, viral capsid antigen.

**TABLE 18** Epstein-Barr virus antibodies and the timing of infections

Possible results				
VCA IgG	VCA IgM	EBNA IgG	EA	Interpretation
–	–	–	–	No previous exposure
+	+	–	±	Acute infection
+	–	+	–	Past infection
+	±	±	±	Recent infection
+	±	+	+	Reactivation

*EA*, Early antigen; *EBNA*, Epstein-Barr virus nuclear antigen; *IgG*, immunoglobulin G; *VCA*, viral capsid antigen.

- Anti-EA (EA-D or EA-R) antibodies appear next or are present with anti-VCA antibodies early in the course of illness. An anti-EA antibody titer greater than 80 in a patient 2 years after acute infectious mononucleosis indicates chronic EBV syndrome.
- As the patient recovers, anti-VCA IgG and anti-EA antibodies decrease, and anti-EBNA antibodies appear. Anti-EBNA antibodies persist for life and reflect a past infection.
- After the patient is well, anti-VCA IgG and anti-EBNA antibodies are always present but at lower ranges. Occasionally, anti-EA antibodies also may be present after the patient recovers.

## 370 Epstein-Barr virus testing

In an acute infection, heterophile antibodies usually appear on Monospot within the first 3 weeks of illness but then decline rapidly within a few weeks. The heterophile antibody, however, fails to develop in about 10% of adults, more frequently in children, and almost uniformly in infants with primary EBV infections. If EBV infection is suspected to have occurred more than a few weeks before testing, the Monospot test result may be negative. Detecting anti-VCA IgG or EBNA will not be helpful because they indicate that an EBV infection has occurred sometime in the patient's life but not necessarily recently. However, detecting anti-VCA IgM would indicate that the syndrome of complaints the patient experienced a few weeks prior was because of EBV.

In immunosuppressed patients (i.e., AIDS, transplantation, or long-term chemotherapy), EBV infection can be much more serious, instigating extranodal lymphoma and posttransplant lymphoproliferative disorders. These patients may have serologic negative tests because of their immunosuppression.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: verify with laboratory
- Obtain serum samples as soon as possible after the onset of the illness.
- Obtain a second blood specimen 14 to 21 days later.

### Abnormal findings

Infectious mononucleosis  
Chronic fatigue syndrome  
Chronic EBV carrier state  
Burkitt lymphoma  
Nasopharyngeal cancer  
Posttransplant lymphoproliferative disease

### notes

**erythrocyte fragility** (Osmotic fragility [OF], Red blood cell fragility)**Type of test** Blood**Normal findings**

Hemolysis begins at 0.5% NaCl.

Hemolysis complete at 0.3% NaCl.

**Test explanation and related physiology**

Red blood cells (RBCs) are bound by a membrane that allows water to pass through while generally restricting the solutes. This process, called *osmosis*, causes RBCs to absorb water when in a hypotonic medium. This results in swelling and, ultimately, hemolysis when the cell bursts. The osmotic fragility (OF) test uses this fact to determine the concentration of solute inside the cell by subjecting it to salt solutions of different concentrations. The ability of the normal RBC to withstand hypotonicity results from its biconcave shape, which allows the cell to increase its volume by 70% before the surface membrane is stretched. When this limit is reached, lysis occurs. When intravascular hemolysis is identified, OF is used to determine whether the RBCs have increased fragility (tend to burst open when exposed to a higher concentrated NaCl solution) or decreased fragility (tend to burst open in lower concentrated, and thus more hypotonic, NaCl solution).

This test is performed to detect hereditary spherocytosis and thalassemia when intravascular hemolysis is identified. Round cells (spherocytes) have increased OF compared with normal indented RBCs. In hereditary spherocytosis, there is abnormal morphology because of a lack of spectrin, a key RBC cytoskeletal membrane protein. This produces membrane instability, which forces the cell to the smallest volume—that of a sphere. This common disorder is associated with intravascular hemolysis. This is shown by increased osmotic fragility. Thalassemia, on the other hand, is associated with thinner leptocytes whose OF is decreased.

**Interfering factors**

- Acute hemolysis. The osmotically labile cells are already hemolyzed and therefore are not found in the blood specimen.
- ✦ Dapsone can increase OF.

## 372 erythrocyte fragility

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: green

### Abnormal findings

▲ Increased erythrocyte	▼ Decreased erythrocyte fragility
Acquired hemolytic anemia	Thalassemia
Hereditary spherocytosis	Hemoglobinopathies (C and S disease)
Hemolytic disease of the newborn	Iron deficiency anemia
Pyruvate kinase deficiency	Reticulocytosis
Malaria	

### notes



## erythrocyte sedimentation rate (ESR, Sed rate test)

**Type of test** Blood

### Normal findings

#### Westergren method

Male:  $\leq 15$  mm/hr

Female:  $\leq 20$  mm/hr

Child:  $\leq 10$  mm/hr

Newborn: 0-2 mm/hr

### Test explanation and related physiology

ESR is a measurement of the rate with which the RBCs settle in saline or plasma over a specified time period. It is nonspecific and therefore not diagnostic for any particular organ disease or injury. Because acute and chronic infection, inflammation (collagen-vascular diseases), advanced neoplasm, and tissue necrosis or infarction increase the protein (mainly fibrinogen) content of plasma, RBCs have a tendency to stack up on one another, increasing their weight and causing them to descend faster. Therefore in these diseases, the ESR is increased. ESR is considered an *acute-phase* or a *reactant* protein (i.e., it occurs as a reaction to acute illnesses).

The test can be used to detect occult disease. Many physicians use the ESR test in this way for routine patient evaluation for vague symptoms. Other physicians regard this test as so nonspecific that it is useless as a routine study. The ESR test occasionally can be helpful in differentiating disease entities or complaints.

The ESR is a fairly reliable indicator of the course of disease and can be used to monitor disease therapy, especially for inflammatory autoimmune diseases (e.g., temporal arteritis or polymyalgia rheumatica). In general, as the disease worsens, the ESR increases; as the disease improves, the ESR decreases. If the results of the ESR are equivocal or inconsistent with clinical impressions, the C-reactive protein test (see p. 295) is often performed.

### Interfering factors

- Artificially low results can occur when the collected specimen is allowed to stand longer than 3 hours before the testing.
- Pregnancy (second and third trimesters) can cause elevations.
- Menstruation can cause elevated levels.
- Polycythemia is associated with decreased ESR.

### 374 erythrocyte sedimentation rate

- ❖ Drugs that may cause *increased* ESR levels include dextran, methyldopa, oral contraceptives, penicillamine, procainamide, theophylline, and vitamin A.
- ❖ Drugs that may cause *decreased* levels include aspirin, cortisone, and quinine.

#### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender (verify with laboratory)

#### Abnormal findings

##### ▲ Increased levels

Chronic renal failure  
Malignant diseases  
Bacterial infection  
Inflammatory diseases  
Necrotic tissue diseases  
Hyperfibrinogenemia  
Macroglobulinemia  
Severe anemias (e.g., iron deficiency or B<sub>12</sub> deficiency)

##### ▼ Decreased levels

Sickle cell anemia  
Spherocytosis  
Hypofibrinogenemia  
Polycythemia vera

notes

## erythropoietin (EPO)

**Type of test** Blood

### Normal findings

5-35 IU/L

### Test explanation and related physiology

EPO is a hormone produced by the kidney. In response to decreased oxygen, the production of EPO is increased. EPO stimulates the bone marrow to increase RBC production. This improves oxygenation in the kidney, and the stimulus for EPO is reduced. This feedback mechanism is very sensitive to minimal persistent changes in oxygen levels. In patients with normal renal function, EPO levels are inversely proportional to the hemoglobin concentration.

As a hormone, EPO is often administered to patients who experience anemia as a result of chemotherapy. Occasionally athletes, in order to improve oxygen-carrying capacity and thereby improve performance, abuse this hormone.

EPO testing is performed to assist in the differential diagnosis of patients with anemia or polycythemia. EPO is elevated in patients who have low hemoglobin due to failure of marrow production or with increased RBC destruction (iron-deficiency or hemolytic anemia, respectively). However, patients with anemia caused by renal diseases (or bilateral nephrectomy) do not have elevated EPO levels. The renal cells are damaged by disease. EPO levels fall and these patients develop anemia.

Patients who have polycythemia as an appropriate response to chronic hypoxemia have elevated EPO levels. Yet patients who have malignant marrow causes of polycythemia vera may have reduced EPO levels. Some renal cell or adrenal carcinomas can produce elevated EPO levels.

### Interfering factors

- Pregnancy is associated with elevated EPO levels.
- The use of transfused blood decreases EPO levels.
- ✦ Drugs that may *increase* EPO levels include adrenocorticotropic hormone, oral contraceptives, and steroids.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or gel separator

### Abnormal findings

#### ▲ Increased levels

- Iron-deficiency anemia
- Megaloblastic anemia
- Hemolytic anemia
- Myelodysplastic syndrome
- Chemotherapy
- AIDS
- Pheochromocytoma
- Renal cell carcinoma
- Adrenal carcinoma

#### ▼ Decreased levels

- Polycythemia vera
- Renal diseases
- Renal failure

### notes

## esophageal function studies (Esophageal manometry, Esophageal motility studies)

**Type of test** Manometric

### Normal findings

Lower esophageal sphincter pressure: 10-20 mm Hg

Swallowing pattern: normal peristaltic waves

Acid reflux: negative

Acid clearing: < 10 swallows

Bernstein test: negative

### Test explanation and related physiology

Esophageal function studies include:

- Determination of the *lower esophageal sphincter (LES)* pressure (manometry)
- Graphic recording of esophageal swallowing waves or *swallowing pattern* (manometry)
- Detection of reflux of gastric acid back into the esophagus (acid reflux)
- Detection of the ability of the esophagus to clear acid (acid clearing)
- An attempt to reproduce symptoms of heartburn (Bernstein test)

### Manometry studies

Two manometry studies are used in assessing esophageal function: measurement of LES pressure and graphic recording of swallowing waves (motility). The LES is a sphincter muscle that acts as a valve to prevent reflux of gastric acid into the esophagus. Free reflux of gastric acid occurs when sphincter pressures are low. An example of such a disorder in adults is gastroesophageal reflux; in children, it is called achalasia (incompetent or relaxed LES).

With increased sphincter pressure, as found in patients with achalasia (failure of the LES to relax normally with swallowing) and diffuse esophageal spasms, food cannot pass from the esophagus into the stomach. Increased LES pressures are noted on manometry.

### Acid reflux with pH probe

Acid reflux is the primary component of gastroesophageal reflux. Patients with an incompetent LES will regurgitate gastric acid into the esophagus. This then causes a drop in the esophageal pH during *esophageal pH monitoring*.

## 378 esophageal function studies

For this test, a thin, small tube with a device on the tip that senses acid is gently passed through a nostril into the esophagus as the patient swallows. The catheter protruding from the nose is connected to a recorder that registers acid reflux and marks meals, sleep, and symptoms. The recorder is worn on a belt or over the shoulder for 24 to 48 hours.

A newer, *wireless pH probe* uses a capsule and makes pH monitoring easier. The capsule is introduced into the esophagus on a catheter through the nose or mouth and attached to the lining of the esophagus. The catheter is then detached from the capsule and removed. No catheter protrudes from the nose. Patients can eat and drink normally as well as engage in their usual activities while having their pH levels tested. Within days, the capsule spontaneously sloughs off the wall of the esophagus and passes into the stool. It is not reusable. After the study is completed, the patient returns the recorder, and the data are downloaded to a computer for analysis.

The use of esophageal electrical impedance studies (EEI) is becoming increasingly performed along with pH monitoring. With the use of a multichannel intraluminal impedance (MII) probe, reflux into the esophagus (for as little as 15 seconds) can be detected. Measuring impedance at multiple sites (multichannel—MMI) allows for determination of the direction of bolus movement. Bolus entries progressing from proximal to distal indicate antegrade bolus movement, while bolus entries progressing from distal to proximal indicate retrograde bolus movement.

### Acid clearing

Patients with normal esophageal function can completely clear hydrochloric acid from the esophagus in less than 10 swallows. Patients with decreased esophageal motility (frequently caused by severe esophagitis) require a greater number of swallows to clear the acid.

### Bernstein test (acid perfusion)

The Bernstein test is simply an attempt to reproduce the symptoms of gastroesophageal reflux. If the patient has pain with the instillation of hydrochloric acid into the esophagus, the test result is positive and proves that the patient's symptoms are caused by reflux esophagitis. If the patient has no discomfort, a cause other than esophageal reflux must be sought to explain the patient's symptoms.

## Contraindications

- Patients who cannot cooperate
- Patients who are medically unstable

## Potential complications

- Aspiration of gastric contents

## Interfering factors

- Eating shortly before the test may affect results.
- ✚ Drugs such as sedatives can alter test results.


## Procedure and patient care

### Before

- PT** Explain the procedure to the patient. Inform the patient that some medications may be stopped up to 7 days before testing. Verify all medications with the physician.
- PT** Instruct the patient not to eat or drink anything for at least 8 hours before the test.
- PT** Allay any fears and allow the patient to verbalize concerns. Be sensitive to the patient's fears about choking.

### During

- Note the following procedural steps:
  1. Esophageal studies are usually performed in the endoscopy laboratory.
  2. The fasting, unседated patient is asked to swallow two or three very tiny tubes. The tubes are equipped so that pressure measurements can be taken at 5-cm intervals ([Figure 20](#)).
  3. The outer ends of the tubes are attached to a pressure transducer.
  4. All tubes are passed into the stomach; then three tubes are slowly pulled back into the esophagus. A rapid and extreme increase in the pressure readings indicates the high pressure zone of the LES.
  5. The LES pressure is recorded.
  6. With all tubes in the esophagus, the patient is asked to swallow. Motility wave patterns are recorded.
  7. The pH indicator probe is placed in the esophagus. As described earlier, this can be done wirelessly over a 24- to 48-hour period.
  8. Acid clearing can be determined by asking the patient to swallow after an acid solution is instilled into the esophagus. More than 10 swallows to clear the acid indicates decreased esophageal motility.

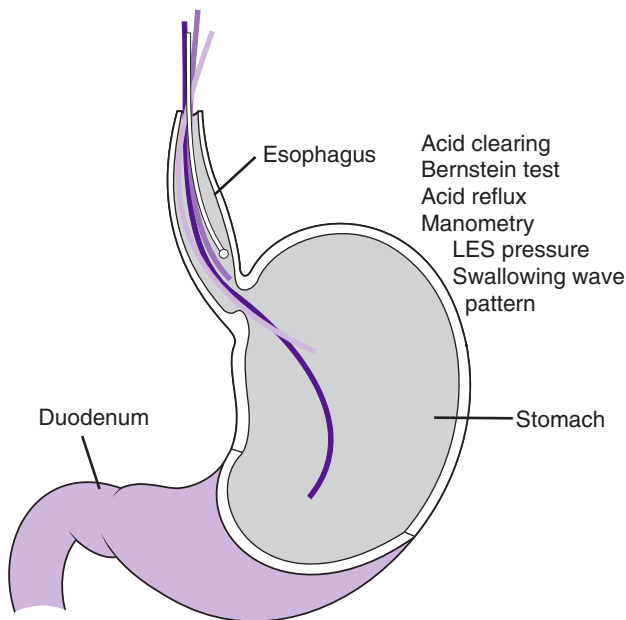

**FIG. 20** Esophageal function studies demonstrating placement of manometry tubes and a pH probe in the esophagus. *LES*, Lower esophageal sphincter.

9. For the Bernstein test, 0.1-N hydrochloric acid and saline solution are alternately instilled into the esophagus. The patient is not told which solution is being infused. If the patient volunteers symptoms of discomfort while the acid is running, the test result is considered positive. If no discomfort is recognized, the test result is negative.

- Note that these tests are performed by an esophageal technician in approximately 30 minutes.
- PT** Inform the patient that the test results are interpreted by a physician and will be available in a few hours.
- PT** Tell the patient that, except for some initial gagging when swallowing the tubes, these tests are not uncomfortable.

#### After

- PT** Inform the patient that it is not unusual to have a mild sore throat after placement of the tubes.



## Abnormal findings

Presbyesophagus  
Diffuse esophageal spasm  
Chalasia  
Achalasia  
Gastroesophageal reflux  
Reflux esophagitis

notes

**esophagogastroduodenoscopy** (EGD, Upper gastrointestinal [UGI] endoscopy, Gastroscopy)**Type of test** Endoscopy**Normal findings**

Normal esophagus, stomach, and duodenum

**Test explanation and related physiology**

Endoscopy enables direct visualization of the upper GI tract by means of a long, flexible, fiberoptic-lighted scope. The esophagus, stomach, and duodenum are examined for tumors, varices, mucosal inflammations, hiatal hernias, polyps, ulcers, and obstructions. It is used to evaluate patients with dysphagia, weight loss, early satiety, upper abdominal pain, ulcer symptoms, or dyspepsia. It is also used to detect esophageal varices in cirrhotic patients. Suspicious barium swallow or upper GI x-ray findings can be corroborated by EGD.

With upper GI endoscopy, one can also visualize and perform a biopsy of tissue in most of the small intestinal tract. This procedure is referred to as *enteroscopy*. Abnormalities of the small intestine such as arteriovenous (AV) malformations, tumors, enteropathies (e.g., celiac disease), and ulcerations can be diagnosed with enteroscopy.

*Capsule endoscopy* (or *wireless capsule endoscopy*) uses a capsule containing a miniature camera that records images of the entire digestive tract, particularly the small intestine. The most common reason for doing capsule endoscopy is to search for a cause of bleeding from the small intestine. It may also be useful for detecting polyps, inflammatory bowel disease (Crohn disease), ulcers, and tumors of the small intestine.

An experienced endoscopist often can control active GI bleeding by electrocoagulation, laser coagulation, or the injection of sclerosing agents (e.g., alcohol for esophageal varices). Also, with the endoscope, benign and malignant strictures can be dilated to reestablish patency of the upper GI tract. Biliary, esophageal, and duodenal stents and percutaneous gastrostomy can be placed with the use of EGD.

**Contraindications**

- Patients with severe upper GI bleeding. The viewing lens will become covered with blood clots.
- Patients with esophageal diverticula. The scope can easily fall into the diverticulum and perforate the wall of the esophagus.

- Patients with suspected perforation. The perforation can be worsened by the insufflation of pressurized air into the GI tract.
- Patients who have had recent GI surgery. The anastomosis may not be able to withstand the pressure of the required air insufflation.

### Potential complications

- Perforation of the esophagus, stomach, or duodenum
- Bleeding from a biopsy site or scope trauma
- Pulmonary aspiration of gastric contents
- Oversedation from medication administered during the test
- Hypotension induced by the sedative medication

### Interfering factors

- Food in the stomach
- Excessive GI bleeding

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- Obtain informed consent if required by the institution.
- PT** Instruct the patient to abstain from eating for 8 to 10 hours before the test.
- If capsule endoscopy is to be performed, apply the recording mechanism to the patient's abdomen and waist.
- PT** Reassure the patient that this test is not painful but may be mildly uncomfortable. Tell the patient that the throat will be anesthetized with a spray to depress the gag reflex.
- PT** Encourage the patient to verbalize fears. Provide support.
- Remove the patient's dentures and eyewear before testing.
- PT** Remind the patient that he or she will not be able to speak during the test but that respiration will not be affected.
- PT** Instruct the patient not to bite down on the endoscope.
- PT** Instruct the patient as to appropriate oral hygiene because the tube will be passed through the mouth.

#### During

- Note the following procedural steps:
  1. The patient is placed on the endoscopy table in the left lateral decubitus position.
  2. The throat is topically anesthetized with viscous lidocaine or another anesthetic spray.

## 384 esophagogastroduodenoscopy

3. The patient is usually sedated. This minimizes anxiety and allows the patient to experience a light sleep.
  4. The endoscope is gently passed through the mouth and finally into the esophagus.
  5. Air is insufflated to distend the upper GI tract for adequate visualization.
  6. The esophagus, stomach, and duodenum are evaluated.
  7. During enteroscopy, the upper small bowel is visualized and a biopsy is performed if needed.
  8. Biopsy or any endoscopic intervention is performed with direct visualization.
  9. On completion of direct inspection and surgery, the excess air and GI secretions are aspirated through the scope.
- Note that the test is performed in the endoscopy laboratory by a physician trained in GI endoscopy and takes approximately 20 to 30 minutes.
- PT** Instruct the patient that eating is allowed about 2 to 4 hours after swallowing the capsule for capsule endoscopy.
- PT** For capsule endoscopy, instruct the patient to return to the physician's office in 6 to 10 hours to return the recording device and have all recording wires removed.
- PT** For capsule endoscopy, instruct the patient that he or she can continue regular activities throughout the examination and will not feel any sensations resulting from the capsule's passage.

### After

- PT** Inform the patient that he or she may have hoarseness or a sore throat after the test. A soothing mouthwash may help.
- Withhold any fluids until the patient is alert and the swallowing reflex returns to normal, usually in 2 to 4 hours.
  - Observe the patient's vital signs. Evaluate the patient for bleeding, fever, abdominal pain, dyspnea, or dysphagia.
  - Observe safety precautions until the effects of the sedatives have worn off.
- PT** For capsule endoscopy, instruct the patient that there is no need to retrieve the capsule or camera from the stool.
- PT** Inform the patient that it is normal to have some bloating, belching, and flatulence after the procedure.
- PT** Inform the patient that the sedation may cause some retrograde and antegrade amnesia for a few hours.

**Abnormal findings**

Tumor (benign or malignant) of the esophagus, stomach, or duodenum

Esophageal diverticula

Hiatal hernia

Esophagitis, gastritis, duodenitis

Gastroesophageal varices

Peptic ulcer

Peptic stricture and subsequent scarring

Extrinsic compression by a cyst or tumor outside the upper GI tract

Source of upper GI bleeding

*Helicobacter pylori* infection

**notes**

## estrogen fractions (Estriol excretion, Estradiol, Estrone)

**Type of test** Urine (24-hour); blood

### Normal findings

	Serum	Urine mcg/24 hours
<i>Estradiol</i>		
Child < 10 years old	< 15 pg/mL	0-6
Adult male	10-50 pg/mL	0-6
Adult female		
Follicular phase	20-350 pg/mL	0-13
Midcycle peak	150-750 pg/mL	4-14
Luteal phase	30-450 pg/mL	4-10
Postmenopause	≤ 20 pg/mL	0-4
<i>Estriol</i>		
Male or child < 10 years old	N/A	1-11
Female, adult		
Follicular phase	N/A	0-14
Ovulatory phase	N/A	13-54
Luteal phase	N/A	8-60
Postmenopausal	N/A	0-11
Female, pregnant		
1st trimester	< 38 ng/mL	0-800
2nd trimester	38-140 ng/mL	800-12,000
3rd trimester	31-460 ng/mL	5000-12,000
<i>Total estrogen</i>		
Male or child < 10 years old	N/A	4-25
Female, nonpregnant	N/A	4-60
Female, pregnant		
1st trimester	N/A	0-800
2nd trimester	N/A	800-5000
3rd trimester	N/A	5000-50,000

### Possible critical values

Values 40% lower than the average of two previous values demand immediate evaluation of fetal well-being during pregnancy.

## Test explanation and related physiology

There are three major estrogens. Estradiol ( $E_2$ ), the most potent estrogen, is produced predominantly in the ovary. Follicle-stimulating hormone (FSH) and luteinizing hormone (LH) stimulate the ovary to produce  $E_2$ , which peaks during the ovulatory phase of the menstrual cycle. This hormone is measured most often to evaluate menstrual and fertility problems, menopausal status, sexual maturity, gynecomastia, feminization syndromes, or as a tumor marker for patients with certain ovarian tumors.

Estrone ( $E_1$ ) is also secreted by the ovary, but most is converted from androstenedione in peripheral tissues.  $E_1$  is the major circulating estrogen after menopause.

Estriol ( $E_3$ ) is the major estrogen in the pregnant female. Serial urine and blood studies for  $E_3$  excretion provide objective means of assessing placental function and fetal normality in high-risk pregnancies. Unfortunately, only severe placental distress will decrease urinary  $E_3$  sufficiently to reliably predict fetoplacental stress. Furthermore, plasma and urinary  $E_3$  levels are normally associated with significant daily variation, which may confuse serial results. Most clinicians use nonstress fetal monitoring (see p. 415) to indicate fetoplacental health.

$E_3$  excretion studies can be done using urine tests or blood studies. A serially increasing estriol/creatinine ratio is a favorable sign in pregnancy. The advantage of the plasma  $E_3$  determination is that it is more easily obtained than a urine specimen and less affected by medications.  $E_3$  is one of the components of the “quad screen” that is obtained in the second trimester of pregnancy to screen for Down syndrome.

## Interfering factors

- Recent administration of radioisotopes may alter test results.
- Glycosuria and urinary tract infections can increase urine  $E_3$  levels.
- ✎ Drugs that may *increase* levels include adrenocorticosteroids, ampicillin, estrogen-containing drugs, phenothiazines, and tetracyclines.
- ✎ Drugs that may *decrease* levels include clomiphene.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- See inside front cover for Routine Urine Testing for spot urine.

### Abnormal findings

#### ▲ Increased levels

- Feminization syndromes
- Precocious puberty
- Ovarian tumor
- Testicular tumor
- Adrenal tumor
- Normal pregnancy
- Hepatic cirrhosis
- Hepatic necrosis
- Hyperthyroidism

#### ▼ Decreased levels

- Failing pregnancy
- Turner syndrome
- Hypopituitarism
- Primary and secondary hypogonadism
- Stein-Leventhal syndrome
- Menopause
- Anorexia nervosa

#### notes



**estrogen receptor assay** (ER assay, ERA, Estradiol receptor)**Type of test** Microscopic examination**Normal findings**

Immunohistochemistry

Negative: &lt; 5% of the cells stain for receptors

Positive: &gt; 5% of the cells stain for receptors

Reverse-transcriptase polymerase chain reaction (RT-PCR)

Negative: &lt; 6.5 units

Positive: &gt; 6.5 units

**Test explanation and related physiology**

The ER assay is useful in determining the prognosis and treatment of breast cancer. The assay is used to determine whether a tumor is likely to respond to endocrine therapy (e.g., tamoxifen, estrogens, aromatase inhibitors, oophorectomy). Tumors with a positive ER assay are more than twice as likely to respond to endocrine therapy than ER-negative tumors.

Specimens are obtained from surgical specimens by a pathologist. Results are usually available in about 1 week.

**Interfering factors**

- Delay in tissue fixation or too long in fixative solution may cause deterioration of receptor proteins and produce lower values.

**Procedure and patient care****Before****PT** Explain the procedure to the patient.**PT** Instruct the patient to discontinue hormones before breast biopsy is performed.

- Before biopsy, a gynecologic history is obtained, including menopausal status and exogenous hormone use.

**During**

- The surgeon obtains tumor tissue.
- This tissue should be placed on ice or in formalin.
- Part of the tissue is used for routine histology. A portion of the paraffin block or a slide containing cancer is used for immunohistochemistry (IHC) staining.

### **390 estrogen receptor assay**

#### **After**

- Results are usually available in 1 week.

#### **Abnormal findings**

Not applicable

**notes**

**ethanol** (Ethyl alcohol, Blood alcohol, Blood EtOH)

**Type of test** Blood; urine; gastric; breath

**Normal findings**

Blood: 0-50 mg/dL, or 0%-0.05%

**Possible critical values**

Blood: > 300 mg/dL

**Test explanation and related physiology**

This test is usually performed to evaluate alcohol-impaired driving or alcohol overdose. Proper collection, handling, and storage of blood alcohol are important for medicolegal cases involving sobriety. Legal testing must be done by specially trained people and must have a strict chain of custody (a paper trail that records sample movement and handling).

Samples tested for legal purposes may include blood, breath, urine, or saliva. Blood is the specimen of choice. Blood is taken from a peripheral vein in living patients and from the aorta in cadavers. Results are given as mg/dL, g/100 mL, or a percentage. Each represents the same amount of alcohol. Blood alcohol concentrations (BACs) greater than 80 mg/dL (0.08%) may cause flushing, slowing of reflexes, and impaired visual activity. Depression of the central nervous system occurs with BACs over 0.1%, and fatalities are reported with levels greater than 0.4%. Persons with BACs less than 0.05% are not considered under the influence of alcohol. Levels greater than 0.05% are considered in most states to be illegal for the operation of motor vehicles and as definite evidence of intoxication.

For legal purposes, when outside of a laboratory or hospital, taking a blood sample for later analysis in the laboratory is not practical or efficient. Breath testing is the most common test performed on automobile drivers. It uses the tail end sample of breath from deep in the lungs and uses a conversion factor to estimate the amount of alcohol in the blood. Blood alcohol testing may be ordered to confirm or refute findings. Alcohol that a person drinks shows up in the breath because it gets absorbed from the intestinal tract into the bloodstream. The alcohol is not metabolized on first pass through the liver. As the blood goes through the lungs, some of the volatile alcohol moves across the alveolar membranes and is exhaled.

Urine testing may also be performed as an alternative to blood. Usually, a patient collects and discards a urine sample and then collects a second sample 20 to 30 minutes later. Saliva alcohol testing is not as widely used but may be used as an alternate screening test. Alcohol stays in the saliva for 6 to 12 hours. Finally, hair testing is used but represents a more chronic use of alcohol.

### **Interfering factors**

- Elevated blood ketones (as with diabetic ketoacidosis) can cause false elevations of blood and breath test results.
  - Bacteria in the urine of diabetic patients with glucosuria can metabolize the glucose to alcohol.
  - Alcohols other than ethanol (e.g., isopropyl [rubbing alcohol] or methanol [grain alcohol]) will cause positive results.
- PT** The use of alcohol-based mouthwash or cough syrup may cause false positives on a breath test.

### **Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: gray or red
- Follow the institution's chain-of-custody protocol.
- Follow the agency's protocol regarding specimen collection.
- Use a povidone-iodine wipe instead of an alcohol wipe for cleansing the venipuncture site.
- If a gastric or urine specimen is indicated, approximately 20 to 50 mL of fluid is necessary.
- Breath analyzers are taken at the end of expiration after a deep inspiration.
- The exact time of specimen collection should be indicated.

### **Abnormal findings**

Alcohol intoxication or overdose

notes

## evoked potential studies (EP studies, Evoked brain potentials, Evoked responses)

**Type of test** Electrodiagnostic

### Normal findings

No neural conduction delay

### Test explanation and related physiology

EP studies are indicated for patients who are suspected of having a sensory deficit but are unable to or cannot reliably indicate recognition of a stimulus. These may include infants, comatose patients, or patients with an inability to communicate. These tests are also used to evaluate specific areas of the cortex that receive incoming stimuli from the eyes, ears, and lower or upper extremity sensory nerves. They are used to monitor natural progression or treatment of deteriorating neurologic diseases.

EP studies focus on changes and responses in brain waves that are evoked from stimulation of a sensory pathway. The EP study measures minute voltage changes produced in response to a specific stimulus such as a light pattern, a click, or a shock. In contrast to the EEG, which records signals that reach amplitudes of up to 50 to 100 millivolts (mV), EP signals are usually less than 5 mV.

Evoked potential studies measure and assess the entire sensory pathway from the peripheral sensory organ all the way to the brain cortex (recognition of the stimulus). Clinical abnormalities are usually detected by an increase in *latency*, which refers to the delay between the stimulus and the wave response. Sensory stimuli used for the EP study can be visual, auditory, or somatosensory. The sensory stimulus chosen depends on which sensory system is suspected of being pathologic (e.g., questionable blindness, deafness, or numbness) and the area of the brain in which pathology is suspected (auditory stimuli check the brainstem and temporal lobes of the brain; visual stimuli test the optic nerve, central neural visual pathway, and occipital parts of the brain; and somatosensory stimuli check the peripheral nerves, spinal cord, and parietal lobe of the brain). Increased latency indicates pathology of the sensory organ or the specific neural pathway as described earlier.

*Visual-evoked responses (VERs)* are usually stimulated by a strobe light flash, reversible checkerboard pattern, or retinal stimuli. Ninety percent of patients with multiple sclerosis show abnormal latencies of VERs, a phenomenon attributed to demyelination

## 394 evoked potential studies

of nerve fibers. In addition, patients with other neurologic disorders (e.g., Parkinson disease) show an abnormal latency of VERs. The degree of latency seems to correlate with the severity of the disease. Abnormal results also may be seen in patients with lesions of the optic nerve, optic tract, visual left, and the eye itself. Absence of binocularity, which is a neurologic developmental disorder in infants, can be detected and evaluated by VERs. Eyesight problems or blindness can be detected in infants through VERs or *electroretinography*. This test can also be used during eye surgery to provide a warning of possible damage to the optic nerve. The gross visual acuity of infants can even be checked via VERs.

*Auditory brainstem-evoked potentials (ABEPs)* are usually stimulated by clicking sounds to evaluate the central auditory pathways of the brainstem. Either ear can be evoked to detect lesions in the brainstem that involve the auditory pathway without affecting hearing. One of the most successful applications of ABEPs has been screening newborns and other infants for hearing disorders. Recognition of deafness enables infants to be fitted with corrective devices as soon as possible before learning to speak (to prevent speech pathology). ABEPs also have great therapeutic implications in the early detection of posterior fossa brain tumors.

*Somatosensory-evoked responses (SERs)* are usually stimulated by sensory stimulus to an area of the body. The time is then measured for the current of the stimulus to travel along the nerve to the cortex of the brain. SERs are used to evaluate patients with spinal cord injuries and to monitor spinal cord functioning during spinal surgery. They are also used to monitor treatment of diseases (e.g., multiple sclerosis), to evaluate the location and extent of areas of brain dysfunction after head injury, and to pinpoint tumors at an early stage. These tests can also be used to identify malingering or hysterical numbness (e.g., that latency is normal in these patients despite the fact that patients indicate numbness).

One of the main benefits of EPs is their objectivity because voluntary patient response is not needed. This objectivity makes EPs useful with nonverbal and uncooperative patients. It permits the distinction of organic from psychogenic problems. This is invaluable in settling lawsuits concerning workers' compensation insurance.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to shampoo his or her hair before the test.
- PT** Tell the patient that no fasting or sedation is required.

## During

- Note that the position of the electrode depends on the type of EP study to be done:
    1. *VERs* are measured using electrodes placed on the scalp along the vertex and cortex lobes. Stimulation occurs by using a strobe light, checkerboard pattern, or retinal stimuli.
    2. *ABEPs* are stimulated with clicking noises or tone bursts delivered via earphones. The responses are detected by scalp electrodes placed along the vertex and on each earlobe.
    3. *SERs* are stimulated using electrical stimuli applied to nerves at the wrist (medial nerve) or knee (peroneal nerve). The response is detected by electrodes placed over the sensory cortex of the opposite hemisphere on the scalp.
  - Note that this study is performed by a physician or technician in less than 30 minutes.
- PT** Tell the patient that little or no discomfort is associated with this study.

## After

- If gel was used for the adherence of the electrodes, remove it.

## Abnormal findings

### Prolonged latency for VER

Parkinson disease

Demyelinating diseases (e.g., multiple sclerosis)

Optic nerve damage

Ocular disease or injury

Blindness

Optic tract disease

Occipital lobe tumor or CVA

Absence of binocularity

Visual field defects

Occipital lobe tumor or CVA

Absence of binocularity

Visual field defects

### Prolonged latency for ABEP

Demyelinating diseases (e.g., multiple sclerosis)

Tumor (e.g., acoustic neuroma)

Cerebrovascular accident

Auditory nerve damage

Deafness

## **396 evoked potential studies**

### **Abnormal latency for SER**

Spinal cord injury

Cervical disc disease

Spinal cord demyelinating diseases

Peripheral nerve injury, transection, or disease

Parietal cortical tumor or CVA

### **notes**



## factor V Leiden (FVL, Leiden factor V, Mutation analysis)

**Type of test** Blood

### Normal findings

Negative for factor V Leiden

### Test explanation and related physiology

Factor V is an important factor in reaction 4 (common pathway) of normal hemostasis (see p. 251). The term *factor V Leiden* (FVL) refers to an inherited abnormal form of the gene for factor V. The endogenous anticoagulant protein C (see p. 751) is normally able to break down factor V at one of these cleavage sites. However, protein C cannot inactivate this same cleavage site on FVL. Therefore FVL is inactivated at a rate approximately 10 times slower than that of normal factor V and persists longer in the circulation. This results in increased thrombin generation and a mild, hypercoagulable state reflected by elevated levels of prothrombin fragment F1 + 2 and other activated coagulation markers. This test is used to diagnose FVL-associated thrombophilia.

Individuals *heterozygous* for the FVL mutation have a slightly increased risk for venous thrombosis. *Homozygous* individuals have a much greater thrombotic risk (e.g., deep vein thrombosis [DVT], arterial thrombosis, or pulmonary embolism).

Individuals who are candidates for FVL testing include patients who

- experienced a thrombotic event without any predisposing factors
- have a strong family history of thrombotic events
- experienced a thrombotic event before 30 years of age
- experienced DVT during pregnancy or while taking birth control pills
- had venous thrombosis at unusual sites (e.g., cerebral, mesenteric, portal, and hepatic veins)
- experienced an arterial clot

Testing for FVL is sometimes preceded by a coagulation screening test called the *activated protein C (APC) resistance test*. This test identifies resistance of factor V to APC. If APC resistance is identified, the patient then may choose to undergo mutation testing by DNA analysis of the *F5* gene, which encodes the factor V protein. This testing should be accompanied by professional genetic counseling for the patient and family members.

## 398 factor V Leiden

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: purple
- If the patient is having FVL mutation analysis, anticoagulants will not interfere with testing.
- As an alternative, genetic testing can be done on the patient's cells obtained by a smear of the oral surface of the cheek.
- Remember: if the patient is receiving anticoagulants, the bleeding time will be increased.

### Abnormal findings

Activated protein C resistance

Factor V Leiden genetic mutation (homozygous or heterozygous)

notes

**febrile antibodies** (Febrile agglutinins)**Type of test** Blood**Normal findings**Titers  $\leq$  1:80**Test explanation and related physiology**

Febrile antibodies are used to support the diagnosis and monitoring of infectious diseases (e.g., salmonellosis, rickettsial diseases, brucellosis, and tularemia). Neoplastic diseases, such as leukemias and lymphomas, are also associated with febrile agglutinins. Appropriate antibiotic treatment of the infectious agent is associated with a drop in the titer activity of febrile antibodies. Rickettsial species produce antibodies that agglutinate *Proteus vulgaris* antigens. This test is nonspecific and insensitive. More specific testing for these infective agents provides more sensitive and specific laboratory testing. Temperature regulation is important for the performance of these tests. Under no circumstances should the febrile agglutinin be heated before delivery to the laboratory.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**Abnormal findings****▲ Increased febrile antibodies**

Salmonellosis infection  
Rickettsial disease  
Brucellosis  
Tularemia leukemia  
Lymphoma  
Systemic lupus erythematosus

**notes**

## fecal calprotectin

**Type of test** Stool

### Normal findings

$\leq 50$  mcg/g

Borderline: 50.1-120 mcg/g

Abnormal:  $\geq 120.1$  mcg/g

### Test explanation and related physiology

Calprotectin comprises a majority of soluble protein content in the cytosol of neutrophils. Elevated fecal calprotectin indicates the migration of neutrophils to the intestinal mucosa that occurs during inflammation. This test is used to identify patients with inflammation of the intestines (e.g., celiac disease) and particularly inflammatory bowel diseases (e.g., Crohn's disease and ulcerative colitis).

The test is helpful as an ancillary diagnostic test for inflammatory bowel diseases and is a biomarker for treatment assessment. It is normal in patients with irritable bowel disease. The level of calprotectin in stool correlates significantly with endoscopic colonic inflammation in both ulcerative colitis and Crohn disease and fecal lactoferrin. Colorectal neoplasia and gastrointestinal (GI) infection also increase fecal calprotectin.

### Interfering factors

✚ Drugs that can increase levels include aspirin and nonsteroidal antiinflammatory medications.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

**PT** Tell the patient that no fasting is needed.

#### During

- Collect a fresh random stool specimen. No preservatives are needed.
- The specimen must be frozen within 18 hours of collection.
- Collect separate specimens when multiple tests are ordered. Do not add to previously collected specimens.
- Note that specimens cannot be collected from a diaper.

## Abnormal findings

### ▲ Increased levels

Crohn disease

Ulcerative colitis

Celiac disease

Infectious colitis

Necrotizing enterocolitis

Antineutrophil cytoplasmic antibody

notes

F

**fecal fat** (Fat absorption, Quantitative stool fat determination)**Type of test** Stool**Normal findings**

Fat: 2-6 g/24 hr or 7-21 mmol/day (SI units)

Retention coefficient:  $\geq 95\%$ **Test explanation and related physiology**

This qualitative or quantitative test is performed to confirm the diagnosis of steatorrhea. Steatorrhea occurs when fat content in the stool is high. It is suspected when the patient has large, greasy, and foul-smelling stools. Determining an abnormally high fecal fat content confirms the diagnosis. Short-gut syndrome and any condition that may cause malabsorption (e.g., sprue, Crohn disease, Whipple disease) or maldigestion (e.g., bile duct obstruction, pancreatic duct obstruction secondary to tumor or gallstones) are also associated with increased fecal fat.

The total output of fecal fat can be tested on a random stool specimen but is more accurate when total 24-, 48-, or 72-hour collection is carried out. Abnormal results from a random specimen should be confirmed by submission of a timed collection. Test values for random fecal fat collections are reported in terms of percent fat.

**Interfering factors**

- ✔ Drugs that may alter test results include enemas and laxatives, especially mineral oil.
- ✔ Drugs that may *decrease* levels of fecal fat include barium and psyllium fiber (e.g., Metamucil).

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to abstain from alcohol ingestion for 3 days before testing.
- PT** Give the patient instructions regarding the appropriate diet (a diet diary may be requested by the laboratory):
  - For adults, usually 100 g of fat per day is suggested for 3 days before and during the collection period.
  - Children, and especially infants, cannot ingest 100 g of fat. Therefore a *fat retention coefficient* is determined by measuring the difference between ingested fat and fecal

fat and then expressing that difference (the amount of fat retained) as a percentage of the ingested fat:

$$\frac{\text{Ingested fat} - \text{Fecal fat}}{\text{Ingested fat}} \times 100\% = \text{Fat retention coefficient}$$

- Note that the normal fat retention coefficient is 95% or greater. A low value indicates steatorrhea.
- For qualitative analysis, the stool is microscopically examined for fat globules.
- PT** Instruct the patient to defecate in a dry, clean container.
- Occasionally, a tongue blade is required to transfer the stool to the specimen container.
- PT** Tell the patient not to urinate in the stool container.
- PT** Inform the patient that even diarrheal stools should be collected.
- PT** Instruct the patient that toilet paper should not be placed in the stool container.
- PT** Tell the patient not to take any laxatives or enemas during this test because they will interfere with intestinal motility and alter test results.

### During

- Collect each stool specimen and send immediately to the laboratory during the 24- to 72-hour testing period. Label each specimen and include the time and date of collection.
- If the specimen is collected at home, give the patient a large stool container to keep in the freezer.

### After

- Inform the patient that a normal diet can be resumed.

## Abnormal findings

### ▲ Increased levels

Cystic fibrosis

Malabsorption secondary to sprue, celiac disease, Whipple disease, Crohn disease, or radiation enteritis

Maldigestion secondary to obstruction of the pancreaticobiliary tree (e.g., cancer, stricture, gallstones)

Short-gut syndrome secondary to surgical resection, surgical bypass, or congenital anomaly

### notes

**ferritin****Type of test** Blood**Normal findings**

Male: 12-300 ng/mL or 12-300 mcg/L (SI units)

Female: 10-150 ng/mL or 10-150 mcg/L (SI units)

Children

Newborn: 25-200 ng/mL

≤ 1 month: 200-600 ng/mL

2-5 months: 50-200 ng/mL

6 months-15 years: 7-142 ng/mL

**Test explanation and related physiology**

The serum ferritin study is a good indicator of available iron stores in the body. Ferritin, the major iron storage protein, is normally present in the serum in concentrations directly related to iron storage. In normal patients, 1 ng/mL of serum ferritin corresponds to approximately 8 mg of stored iron.

Decreases in ferritin levels indicate a decrease in iron storage associated with iron deficiency anemia. A ferritin level less than 10 ng/100 mL is diagnostic of iron deficiency anemia. Increased levels are a sign of iron excess, as seen in hemochromatosis, hemosiderosis, iron poisoning, or recent blood transfusions. Increased ferritin is also noted in patients with megaloblastic anemia, hemolytic anemia, and chronic hepatitis. Furthermore, ferritin is factitiously elevated in patients with chronic disease states, such as neoplasm, alcoholism, uremia, collagen diseases, or chronic liver diseases.

A limitation of this study is that ferritin levels also can act as an acute-phase reactant protein and may be elevated in conditions not reflecting iron stores (e.g., acute inflammatory diseases, infections, metastatic cancer, lymphomas). Elevations in ferritin occur 1 to 2 days after onset of the acute illness and peak at 3 to 5 days.

**Interfering factors**

- Recent transfusions and recent ingestion of a meal with high iron content may cause elevated ferritin levels.
- Recent administration of a radionuclide can cause abnormal levels if testing is performed by radioimmunoassay.
- Hemolytic diseases may be associated with an artificially high iron content.



- Disorders of excessive iron storage (e.g., hemochromatosis, hemosiderosis) are associated with high ferritin levels.
  - Iron-deficient menstruating women may have decreased ferritin levels.
  - Acute and chronic inflammatory conditions and Gaucher disease can falsely increase ferritin levels.
- ▼ Iron preparations may *increase* ferritin levels.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Hemochromatosis  
 Hemosiderosis  
 Megaloblastic anemia  
 Hemolytic anemia  
 Alcoholic or inflammatory  
 hepatocellular disease  
 Inflammatory disease  
 Advanced cancers  
 Chronic illnesses (e.g.,  
 leukemias, cirrhosis,  
 chronic hepatitis)  
 Collagen vascular diseases  
 Hemophagocytic syndromes  
 Congenital and acquired  
 sideroblastic anemias

#### ▼ Decreased levels

Severe protein deficiency  
 Iron-deficiency anemia  
 Hemodialysis

notes

## fetal biophysical profile (Biophysical profile [BPP])

**Type of test** Ultrasound; fetal activity study

### Normal findings

Score of 8-10 points (if amniotic fluid volume is adequate)

### Possible critical values

Score of less than 4 points may necessitate immediate delivery.

### Test explanation and related physiology

The BPP is a method of evaluating fetal status during the antepartal period based on five variables originating within the fetus: fetal heart rate (FHR), fetal breathing movement, gross fetal movements, fetal muscle tone, and amniotic fluid volume. FHR reactivity is measured by the nonstress test (NST) (see p. 415); the other four parameters are measured by ultrasound scanning. BPP is often done to assess fetal well-being in the face of a nonreactive NST.

The major premise behind the BPP is that variable assessments of fetal biophysical activity are more reliable than an examination of a single parameter (e.g., FHR). Indications for this test include such factors as postdate pregnancy, maternal hypertension, diabetes mellitus, vaginal bleeding, maternal Rh sensitization, maternal history of stillbirth, and premature rupture of membranes.

The five parameters are briefly described here. Each parameter is scored and contributes either a 2 or a 0 to the score. Therefore 10 is the perfect score, and 0 is the lowest score. Gestational age influences these results. For example, fetal breathing movements are the latest parameter to develop.

- *FHR reactivity.* This is measured and interpreted in the same way as the nonstress test (see p. 415). The FHR is considered reactive when there are movement-associated FHR accelerations of at least 15 beats/min above baseline and 15 seconds in duration over a 20-minute period. A score of 2 is given for reactivity; a score of 0 indicates that the FHR is nonreactive.
- *Fetal breathing movements.* This parameter is assessed based on the assumption that fetal breathing movements indicate fetal well-being and their absence may indicate hypoxemia. To earn a score of 2, the fetus must have at least one episode of fetal breathing lasting at least 60 seconds within a 30-minute observation. Absence of this breathing pattern is scored a 0 on the BPP. It is important to note that fetal breathing

movements increase during the second and third hours after maternal meals and at night. Fetal breathing movements may decrease in such conditions as hypoxemia, hypoglycemia, nicotine use, and alcohol ingestion.

- *Fetal body movements.* Fetal activity is a reflection of neurologic integrity and function. The presence of at least three discrete episodes of fetal movements within a 30-minute observation period is given a score of 2. A score of 0 is given with two or fewer fetal movements in this time period. Fetal activity is greatest 1 to 3 hours after the mother has consumed a meal.
- *Fetal tone.* In the uterus, the fetus is normally in a position of flexion. However, the fetus also stretches, rolls, and moves in the uterus. The arms, legs, trunk, and head may be flexed and extended. A score of 2 is earned when there is at least one episode of active extension with return to flexion. An example of this is the opening and closing of a hand. A score of 0 is given for slow extension with a return to only partial flexion, fetal movement not followed by return to flexion, limbs or spine in extension, and an open fetal hand.
- *Amniotic fluid volume.* Oligohydramnios (too little amniotic fluid) has been associated with fetal anomalies, intrauterine growth retardation, and postterm pregnancy. A score of 2 is given for this parameter when there is at least one pocket of amniotic fluid that measures 1 cm in two perpendicular planes. A score of 0 indicates either that fluid is absent in most areas of the uterine cavity or that the largest pocket measures 1 cm or less in the vertical axis.

A score of 8 or 10 with a normal amount of amniotic fluid indicates a healthy fetus. A score of 8 with oligohydramnios or a score of 4 to 6 is equivocal. An equivocal test result is interpreted as possibly abnormal. A score of 0 or 2 is abnormal and indicates the need for assessment of immediate delivery.

Modifications can be made to the BPP. Some physicians omit the nonstress test if the ultrasound parameters are normal. Some physicians have added placental grading as a sixth parameter.

Another measure of fetal well-being is the *amniotic fluid index (AFI)*. Ultrasound is used to measure the largest collection of amniotic fluid in each of the four quadrants within the uterus. The numbers are added together, and the sum is the AFI. The normal range for the AFI is 8 to 18 cm. The sum is plotted on a graph in which the age of gestation also is taken into account. If the AFI is less than the 2.4 percentile, oligohydramnios is present. If AFI exceeds the 97th percentile, polyhydramnios exists.

## 408 fetal biophysical profile

An abnormally low AFI observed in antepartum testing is associated with an increased risk of intrauterine growth restriction and overall adverse perinatal outcome. The percentile value seems to be a better indicator than an absolute fluid volume. Oligohydramnios is associated with placental failure or fetal renal problems. Polyhydramnios is associated with maternal diabetes or fetal upper GI malformation or obstruction.

Doppler ultrasound evaluations of the placenta and the *umbilical artery velocity* can recognize alterations in umbilical artery flow and direction that may indicate fetal stress or illness.

### Interfering factors

- If no eye movement or respiratory movement is noted, the fetus may be sleeping. Testing is then extended.
- ✚ Central nervous system (CNS) stimulants (e.g., catecholamine) can *increase* BPP activities.
- ✚ Antenatal steroids can *increase* BPP activities.
- ✚ Sedatives and narcotic analgesics can *decrease* BPP activities.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

**PT** Inform the patient that no fasting is required.

#### During

- FHR reactivity is measured and interpreted from a nonstress test (see p. 415).
- Fetal breathing movements, fetal body movements, fetal tone, and amniotic fluid volume are determined by ultrasound imaging (see obstetric ultrasonography, p. 685).

#### After

- If the test results are abnormal or equivocal, support the patient in the next phase of the fetal evaluation process.

### Abnormal findings

Fetal asphyxia

Postterm pregnancy

Congenital anomalies

Fetal stress

Oligohydramnios

Fetal death

Intrauterine growth retardation

#### notes

**fetal contraction stress test** (Contraction stress test [CST],  
Oxytocin challenge test [OCT])

**Type of test** Electrodiagnostic

**Normal findings**

Negative

**Test explanation and related physiology**

The CST, frequently called the *oxytocin challenge test (OCT)*, is a relatively noninvasive test of fetoplacental adequacy used in the assessment of high-risk pregnancy. For this study, a temporary stress in the form of uterine contractions is applied to the fetus after the IV administration of oxytocin. The reaction of the fetus to the contractions is assessed by an external fetal heart monitor. Uterine contractions cause a transient impediment of placental blood flow. If the placental reserve is adequate, the maternal–fetal oxygen transfer is not significantly compromised during the contractions and the FHR remains normal (a *negative* test result). The fetoplacental unit can then be considered adequate for the next 7 days.

If the placental reserve is inadequate, the fetus does not receive enough oxygen during the contraction. This results in intrauterine hypoxia and late deceleration of the FHR. The test result is considered to be *positive* if consistent, persistent, late decelerations of the FHR occur with two or more uterine contractions. False-positive results caused by uterine hyperstimulation can occur in 10% to 30% of patients. Thus positive test results warrant a complete review of other studies (e.g., amniocentesis) before the pregnancy is terminated by delivery.

The test is considered to be *unsatisfactory* if the results cannot be interpreted (e.g., because of hyperstimulation of the uterus, excessive movement of the mother, or deceleration of unknown meaning). In the case of unsatisfactory results, other means of management should be considered.

Although this test can be performed reliably at 32 weeks of gestation, it usually is done after 34 weeks. CST can induce labor, and a fetus at 34 weeks is more likely to survive an unexpectedly induced delivery than a fetus at 32 weeks. Nonstress testing (see p. 415) of the fetus is the preferred test in almost every instance and can be performed more safely at 32 weeks; it can then be followed 2 weeks later by CST if necessary. The CST may be performed weekly until delivery terminates pregnancy.

## 410 fetal contraction stress test

### Contraindications

- Patient pregnant with multiple fetuses because the myometrium is under greater tension and is more likely to be stimulated to premature labor
- Patient with a prematurely ruptured membrane, because labor may be stimulated by the CST
- Patient with placenta previa because vaginal delivery may be induced
- Patient with abruption placentae because the placenta may separate from the uterus as a result of uterine contractions
- Patient with a previous hysterotomy because the strong uterine contractions may cause uterine rupture
- Patient with a previous vertical or classic cesarean section because the strong uterine contractions may cause uterine rupture.
- Patient with pregnancy of less than 32 weeks because early delivery may be induced by the procedure

### Potential complications

- Premature labor

### Interfering factors

- Hypotension may cause false-positive results.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Obtain informed consent for the procedure.

**PT** Teach the patient breathing and relaxation techniques.

- Record the patient's blood pressure and FHR before the test as baseline values.
- If the CST is performed on an elective basis, the patient may be kept nothing by mouth (NPO) in case labor occurs.

#### During

- Note the following procedural steps:
  1. After the patient empties her bladder, place her in a semi-Fowler position and tilted slightly to one side to avoid vena caval compression by the enlarged uterus.
  2. Check her blood pressure every 10 minutes to avoid hypotension, which may cause diminished placental blood flow and a false-positive test result.
  3. Place an external fetal monitor over the patient's abdomen to record the fetal heart tones. Attach an external

tocodynamometer to the abdomen at the fundal region to monitor uterine contractions.

4. Record the output of the fetal heart tones and uterine contractions on a two-channel strip recorder.
  5. Monitor baseline FHR and uterine activity for 20 minutes.
  6. If uterine contractions are detected during this pretest period, withhold oxytocin and monitor the response of the fetal heart tone to spontaneous uterine contractions.
  7. If no spontaneous uterine contractions occur, administer oxytocin (Pitocin) by intravenous (IV) infusion pump.
  8. Increase the rate of oxytocin infusion until the patient is having moderate contractions; then record the FHR pattern.
  9. After the oxytocin infusion is discontinued, continue FHR monitoring for another 30 minutes until the uterine activity has returned to its preoxytocin state. The body metabolizes oxytocin in approximately 20 to 25 minutes.
- Note that the CST is performed safely on an outpatient basis in the labor and delivery unit, where qualified nurses and necessary equipment are available. The test is performed by a nurse with a physician available.
  - Note that the duration of this study is approximately 2 hours.
- PT** Tell the patient that the discomfort associated with the CST may consist of mild labor contractions. Usually, breathing exercises are sufficient to control any discomfort. Administer analgesics if needed.

### After

- Monitor the patient's blood pressure and FHR.
- Discontinue the IV line and assess the site for bleeding.

### Abnormal findings

Fetoplacental inadequacy

### notes

## fetal fibronectin (fFN, Fibronectin)

**Type of test** Fluid analysis

### Normal findings

Negative:  $\leq 0.05$  mcg/mL

### Test explanation and related physiology

Fibronectin may help with implantation of the fertilized egg into the uterine lining. Normally, fibronectin cannot be identified in vaginal secretions after 22 weeks of pregnancy. However, concentrations are very high in the amniotic fluid. If fibronectin is identified in vaginal secretions after 24 weeks, the patient is at high risk for preterm delivery. The use of fFN is limited to symptomatic women with contractions whose membranes are intact and who have cervical dilation of less than 3 cm.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting is required.
  - Determine whether the patient has had a cervical examination or intercourse within the past 24 hours. Results may be inaccurate.

#### During

- Note the following procedural steps:
  1. The patient is placed in the lithotomy position.
  2. A vaginal speculum is inserted to expose the cervix.
  3. Vaginal secretions are collected from the posterior vagina and paracervical area using a swab from a kit.
  4. The container with appropriate medium is labeled with the patient's name, age, and estimated date of confinement.
- PT** Tell the patient that no discomfort, except for insertion of the speculum, is associated with this procedure.
  - Note that this procedure is performed by a physician or other licensed health care provider in several minutes.

#### After

- PT** Tell the patient that results will be available the same or next day.
- PT** Educate the patient about the signs of preterm labor.

### Abnormal findings

High risk for preterm premature delivery



**fetal hemoglobin** (Kleihauer–Betke)

**Type of test** Blood

**Normal findings**

< 1% of red blood cells (RBCs)

**Test explanation and related physiology**

Fetal hemoglobin may be present in the mother's blood because of fetal–maternal hemorrhage (FMH), which causes leakage of fetal cells into the maternal circulation. When large volumes of fetal blood are lost in this way, neonatal outcomes can be serious and potentially fatal. Massive FMH may be the cause of about 1 in every 50 stillbirths.

Leakage of fetal RBCs can begin any time after the middle of the first trimester. It presumably results from a breach in the integrity of the placental circulation. As pregnancy continues, more women will show evidence of fetal RBCs in their circulation; by term, about 50% will have detectable fetal cells. Most of these, however, are the result of very small leaks. In 96% to 98% of pregnancies, the total fetal blood volume lost in this way is 2 mL or less. Small leaks are not implicated in intrauterine death.

Risk factors correlated with the increasing risk of massive FMH include maternal trauma, placental abruption, placental tumors, third-trimester amniocentesis, fetal hydrops, pale fetal organs, antecedent sinusoidal fetal heart tracing, and twinning. Having one or more of these features should be an indication for fetal hemoglobin testing.

The standard method of detecting FMH is the *Kleihauer–Betke test*. The *flow cytometric method* for fetal hemoglobin determination offers several advantages over the traditional Kleihauer–Betke method. This more objective method has been shown to improve sensitivity, precision, and linearity over traditional methods.

FMH becomes of even greater significance when the mother is Rh negative because this is the mechanism through which Rh sensitization could develop if the fetus has paternal Rh-positive blood cells.

This test is often performed on women who have delivered a stillborn baby to see if FMH was a potential cause of fetal death.

**Interfering factors**

- Any maternal condition (e.g., sickle cell disease) that involves persistence of fetal hemoglobin in the mother will cause a false positive.

## 414 fetal hemoglobin

- If the blood is drawn after C-section, a false positive could occur. Vaginal delivery results in higher frequency of detection of FMH.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Provide emotional support in the event this test is performed after a stillborn delivery.

### Abnormal findings

Fetal–maternal hemorrhage

Hereditary persistence of fetal hemoglobin

Intrachorionic thrombi

notes

## fetal nonstress test (Nonstress test [NST], Fetal activity determination)

**Type of test** Electrodiagnostic

### Normal findings

Reactive fetus (heart rate acceleration associated with fetal movement)

### Test explanation and related physiology

The NST is a method of evaluating the viability of a fetus. It documents the function of the placenta in its ability to supply adequate blood to the fetus. The NST can be used to evaluate any high-risk pregnancy in which fetal well-being may be threatened.

The NST is a noninvasive study that monitors acceleration of the FHR in response to fetal movement. This FHR acceleration reflects the integrity of the CNS and fetal well-being. Fetal activity may be spontaneous, induced by uterine contraction, or induced by external manipulation. Oxytocin stimulation is not used. Fetal response is characterized as *reactive* or *nonreactive*. The NST indicates a reactive fetus when, with fetal movement, two or more FHR accelerations are detected, each of which must be at least 15 beats/min for 15 seconds or more within any 10-minute period. The test is 99% reliable in indicating fetal viability and negates the need for the contraction stress test (CST, p. 409). If the test detects a nonreactive fetus (i.e., no FHR acceleration with fetal movement) within 40 minutes, the patient is a candidate for the CST. A 40-minute test period is used because this is the average duration of the sleep-wake cycle of the fetus. The cycle may vary considerably, however.

An NST is routinely performed before the CST to avoid the complications associated with oxytocin administration. No complications are associated with the NST.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
  - Encourage verbalization of the patient's fears. The need for the study usually raises fears in the expectant mother.
- PT** If the patient is hungry, instruct her to eat before the NST has begun. Fetal activity is enhanced with a high maternal serum glucose level.

## 416 fetal nonstress test

### During

- After the patient empties her bladder, place her in Sims position.
- Place an external fetal monitor on the patient's abdomen to record the FHR. The mother can indicate fetal movement by pressing a button on the fetal monitor whenever she feels the fetus move. FHR and fetal movement are concomitantly recorded on a two-channel strip graph.
- Observe the fetal monitor for FHR accelerations associated with fetal movement.
- If the fetus is quiet for 20 minutes, stimulate fetal activity by external methods, such as rubbing or compressing the mother's abdomen, ringing a bell near the abdomen, or placing a pan on the abdomen and hitting the pan.
- Note that a nurse performs the NST in approximately 20 to 40 minutes in the physician's office or a hospital unit.

**PT** Tell the patient that no discomfort is associated with the NST.

### After

**PT** If the results detect a nonreactive fetus, inform the patient that she is a candidate for the CST. Provide appropriate education.

### Abnormal findings

Fetal stress

Fetal death

### notes

## fetal scalp blood pH

**Type of test** Blood

### Normal findings

pH: 7.25-7.35

O<sub>2</sub> saturation: 30%-50%

Po<sub>2</sub>: 18-22 mm Hg

Pco<sub>2</sub>: 40-50 mm Hg

Base excess: 0 to -10 mEq/L

### Test explanation and related physiology

Measurement of fetal scalp blood pH provides valuable information on fetal acid-base status. This test is useful clinically for diagnosing fetal distress. The pH normally ranges from 7.25 to 7.35 during labor; a mild decline within the normal range is noted with contractions and as labor progresses.

Fetal hypoxia causes anaerobic glycolysis, resulting in excess production of lactic acid. This causes an increase in hydrogen ion concentration (acidosis) and a decrease in pH. Acidosis reflects the effect of hypoxia on cellular metabolism. A high correlation exists between low pH levels and low Apgar scores.

*Fetal oxygen saturation monitoring (FSpo<sub>2</sub>)* also is available to assist the monitoring of fetal well-being during delivery. When the FHR becomes significantly abnormal, C-section is often performed because of concern for fetal well-being. However, with FSpo<sub>2</sub>, an accurate measure of fetal oxygen saturation can be determined. After membranes are ruptured, and if the baby is in vertex position with good cervical dilatation, a specialized probe can be placed on the temple or cheek of the fetus for FSpo<sub>2</sub> monitoring. The oxygen saturation is displayed on a monitor screen as a percentage. The normal oxygen saturation for a baby in the womb receiving oxygenated blood from the placenta is usually between 30% and 70%.

### Contraindications

- Patients with premature membrane rupture
- Patients with active cervical infection

### Potential complications

- Continued bleeding from the puncture site
- Hematoma
- Ecchymosis
- Infection

## 418 fetal scalp blood pH

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Obtain informed consent for this procedure.

**PT** Tell the patient that no fasting or sedation is required.

#### During

- Note the following procedural steps for fetal scalp pH:
    1. Amnioscopy is performed with the mother in the lithotomy position.
    2. The cervix is dilated, and the endoscope (amnioscope) is introduced into the cervical canal.
    3. The fetal scalp is cleansed with an antiseptic and dried with a sterile cotton ball.
    4. A small amount of petroleum jelly is applied to the fetal scalp to cause droplets of fetal blood to bead.
    5. After the skin on the scalp is pierced with a small metal blade, beaded droplets of blood are collected in long, heparinized capillary tubes.
    6. The tube is sealed with wax and placed on ice to retard cellular respiration, which can alter the pH.
    7. The physician performing the procedure applies firm pressure to the puncture site to retard bleeding.
    8. Scalp blood sampling can be repeated as necessary.
  - Note that this study is performed by a physician in approximately 10 to 15 minutes.
- PT** Tell the patient that she may be uncomfortable during the cervical dilation.

#### After

**PT** Inform the patient that she may have vaginal discomfort and menstrual-type cramping.

#### After delivery

- Assess the newborn and document the puncture site(s).
- Cleanse the fetal scalp puncture site with an antiseptic solution and apply an antibiotic ointment.

### Abnormal findings

Fetal distress

notes

## fetoscopy

**Type of test** Endoscopy

### Normal findings

No fetal distress

### Test explanation and related physiology

Fetoscopy is an endoscopic procedure that allows direct visualization of the fetus via the insertion of a tiny, telescope-like instrument through the abdominal wall and into the uterine cavity (Figure 21). Direct visualization may lead to the diagnosis of a severe malformation (e.g., a neural tube defect). During the procedure, fetal blood samples to detect congenital blood disorders (e.g., hemophilia, sickle cell anemia) can be drawn from a blood vessel in the umbilical cord for biochemical analysis. Fetal skin biopsies also can be done to detect primary skin disorders.

Fetoscopy is performed at approximately 18 weeks of gestation. At this time, the vessels of the placental surface are of adequate size, and the fetal parts are readily identifiable. A therapeutic abortion would not be as hazardous at this time as it would be if it were done later in the pregnancy. An ultrasound is usually performed the day after the procedure to confirm the adequacy of the amniotic fluid and fetal viability.

### Potential complications


- Spontaneous abortion
- Premature delivery
- Amniotic fluid leak
- Intrauterine fetal death
- Amnionitis

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Obtain informed consent.
- Assess the FHR to serve as a baseline value.
- Administer fentanyl, if ordered, before the test because it crosses the placenta and quiets the fetus. This prevents excessive fetal movement, which would make the procedure more difficult.

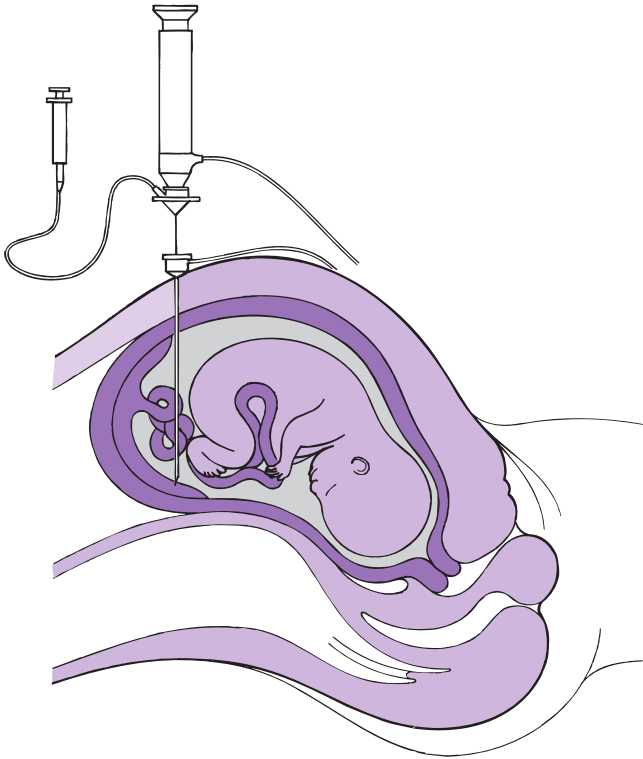

**FIG. 21** Fetoscopy for fetal blood sampling.

### During

- Note the following procedural steps:
  1. The woman is placed in the supine position on an examining table.
  2. The abdominal wall is anesthetized locally.
  3. Ultrasonography is performed to locate the fetus and the placenta.
  4. The endoscope is inserted.
  5. Biopsies and blood samples may be obtained.
- Note that this procedure is performed by a physician in 1 to 2 hours.
- PT** Tell the patient that the only discomfort associated with this study is the injection of the local anesthetic.



### After

- Assess the FHR and compare with the baseline value to detect any side effects related to the procedure.
  - Monitor the mother and fetus carefully for alterations in blood pressure, pulse, uterine activity, and fetal activity; vaginal bleeding; and loss of amniotic fluid.
  - Administer RhoGAM to mothers who are Rh negative unless the fetal blood is found to be Rh negative.
  - Note that a repeat ultrasound is usually performed the day after the procedure to confirm the adequacy of the amniotic fluid and fetal viability.
  - If ordered, administer antibiotics prophylactically after the test to prevent amnionitis.
- PT** Instruct the mother to avoid strenuous activity for 1 to 2 weeks after the procedure.
- PT** Advise the mother to report any pain, bleeding, amniotic fluid loss, or fever.

### Abnormal findings

Developmental defects (e.g., neural tube defects)

Congenital blood disorders (e.g., hemophilia, sickle cell anemia)

Primary skin disorders

### notes

**fibrinogen** (Factor I, Quantitative fibrinogen)**Type of test** Blood**Normal findings**

Adult: 200-400 mg/dL or 2-4 g/L (SI units)

Newborn: 125-300 mg/dL

**Possible critical values**

&lt; 100 mg/dL

**Test explanation and related physiology**

Fibrinogen is essential to the blood-clotting mechanism. It is part of the common pathway in the coagulation system. (See discussion of coagulating factors, p. 248). Fibrinogen, which is produced by the liver, is also an acute-phase protein reactant. It rises sharply during instances of tissue inflammation or necrosis. High levels of fibrinogen have been associated with an increased risk of coronary heart disease, stroke, myocardial infarction, and peripheral arterial disease. This makes fibrinogen an important risk factor for cardiovascular disease.

Fibrinogen is used primarily to aid in the diagnosis of suspected bleeding disorders. This testing is used to detect increased or decreased fibrinogen concentration of acquired or congenital origin. It is also used for monitoring the severity and treatment of disseminated intravascular coagulation and fibrinolysis.

Reduced levels of fibrinogen can be seen in patients with liver disease, malnourished states, and consumptive coagulopathies (e.g., disseminated intravascular coagulation). Large-volume blood transfusions are also associated with low levels of fibrinogen because banked blood does not contain fibrinogen. Reduced levels of fibrinogen will cause a prolonged prothrombin time and partial thromboplastin time.

**Interfering factors**

- Blood transfusions within the past month may affect test results.
- Diets rich in omega-3 and omega-6 fatty acids reduce fibrinogen levels.
- ✎ Drugs that may cause *increased* levels include estrogens and oral contraceptives.
- ✎ Drugs that may cause *decreased* levels include anabolic steroids, androgens, L-asparaginase, phenobarbital, streptokinase, tissue plasminogen activators (e.g., urokinase), and valproic acid.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: blue

## Abnormal findings

### ▲ Increased levels

Acute inflammatory reactions (e.g., rheumatoid arthritis, glomerulonephritis)  
Trauma  
Acute infection (e.g., pneumonia)  
Coronary heart disease  
Cigarette smoking  
Pregnancy  
Cerebrovascular accident  
Myocardial infarction  
Peripheral arterial disease

### ▼ Decreased levels

Liver disease (e.g., hepatitis, cirrhosis)  
Consumptive coagulopathy  
Fibrinolysins  
Congenital afibrinogenemia  
Advanced carcinoma  
Malnutrition  
Large-volume blood transfusion  
Hemophagocytic lymphohistiocytosis

notes

**fluorescein angiography** (FA, Ocular photography)**Type of test** Other**Normal findings**

Normal retinal and choroidal vasculature

**Test explanation and related physiology**

With the use of fluorescein angiography, the patency and integrity of the retinal circulation can be determined. This test involves injection of sodium fluorescein into the systemic circulation followed by timed-interval photographs performed with a fundus camera. The timed images are then reviewed for specific patterns indicative of disease states. This test is performed to diagnose disease affecting the posterior eye, including the retina, choroid, and optic nerve. The test is often repeated at intervals to monitor treatment or disease progression.

Fluorescein enters the ocular circulation from the internal carotid artery and then the ophthalmic artery. Pathologic changes are recognized by the detection of either hyperfluorescence or hypofluorescence. Among the common groups of ophthalmologic disease, fluorescein angiography can detect diabetic retinopathy, vein occlusions, retinal artery occlusions, edema of the optic disc, and tumors. These abnormalities can be treated with a laser to help prevent loss of vision, and treatment results can be monitored using fluorescein angiography.

**Potential complications**

- Allergic reactions. Allergies to fluorescein dye are rare. If they occur, they may cause a skin rash and itching.

**Procedure and patient care****Before****PT** Explain the procedure to the patient.

- Obtain an informed consent.

**PT** Reinforce the need for the patient to remain still during the few seconds after fluorescein injection.

- Obtain an ocular history of cataracts, prior retinal surgery, and other diseases that may inhibit photography.

**PT** Instruct the patient to remove any ocular lenses.**PT** Inform the patient that there are no dietary restrictions.

- Note that pupil dilatation can improve access to the posterior eye. If ordered, administer appropriate mydriatic medications.

Note, however, that these medications are contraindicated for patients with glaucoma because they may dangerously increase ocular pressures.

### During

- Note the following procedural steps:
  1. The patient is positioned in the fundus camera with the chin on the bar.
  2. The patient is told to pick a spot in the far distance and to concentrate on that spot during the examination.
  3. Intravenous access is obtained.
  4. Fluorescein dye is injected with the assistance of an autoinjector.
  5. Photographs are taken by the ophthalmologist at timed intervals.
- The test is performed and interpreted by an ophthalmologist, usually in the office setting. Results are available in less than 30 minutes.

### After

- Remove the intravenous access device and apply pressure to the venipuncture site.
- PT** Inform the patient that fluorescein dye is excreted by the kidneys and to expect very yellow urine for the next 24 hours.

### Abnormal findings

#### ▲ Increased levels

Tumor  
 Detached retina  
 Trauma  
 Inflammation  
 Retinitis pigmentosa  
 Papilledema  
 Diabetic retinopathy

#### ▼ Decreased levels

Diabetes  
 Vascular disease  
 Radiation to the eye  
 Hemorrhage  
 Edema  
 Prior photocoagulation therapy

**folic acid (Folate)****Type of test** Blood**Normal findings**

Serum: 5-25 ng/mL or 11-57 nmol/L (SI units)

RBC folate: 360-1400 nmol/L

**Test explanation and related physiology**

Folic acid (folate), one of the B vitamins, is necessary for normal function of RBCs and white blood cells (WBCs). It is needed for the adequate synthesis of certain purines and pyrimidines, which are precursors for DNA. Blood folate levels require normal absorption by the intestinal tract. The finding of a low serum folate level means that the patient's recent diet has been subnormal in folate content and/or that recent absorption of folate has been subnormal.

Tissue folate is best tested by determining the content of folate in RBCs. A low RBC folate can mean either that there is tissue folate depletion because of folate deficiency, which requires folate therapy, or that the patient has primary vitamin B<sub>12</sub> (see p. 969) deficiency, which blocks the ability of cells to take up folate. In the latter case, the proper therapy would be with vitamin B<sub>12</sub> rather than with folic acid.

Folic acid blood levels are performed to assess folate availability in pregnancy, to evaluate hemolytic disorders, and to detect anemia caused by folic acid deficiency (in which the RBCs are abnormally large, causing a megaloblastic anemia). These RBCs have a shortened life span and impaired oxygen-carrying capacity. If folic acid blood levels are low, RBC folate is measured.

Folate deficiency is present in about 33% of pregnant women; many people with alcoholism; and patients with a variety of malabsorption syndromes, including celiac disease, sprue, Crohn disease, and jejunal/ileal bypass procedure. Folate binds to aluminum hydroxide. Patients with a chronic use of antacids or H<sub>2</sub>-receptor antagonists and with diets marginal in folate may experience low folate levels.

Elevated serum levels of folic acid may be seen in patients with pernicious anemia because vitamin B<sub>12</sub> is needed to allow incorporation of folate into tissue cells. The folic acid tests are often done in conjunction with tests for vitamin B<sub>12</sub> levels.

## Interfering factors

- A falsely normal result may occur in a folate-deficient patient who has received a blood transfusion.
- ✚ Drugs that may cause *decreased* folic acid levels include alcohol, aminopterin, aminosalicic acid, antimalarials, chloramphenicol, erythromycin, estrogens, methotrexate, oral contraceptives, penicillin derivatives, phenobarbital, phenytoins, and tetracyclines.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: verify with laboratory
  - Blood tube commonly used: red
  - Some laboratories prefer an 8-hour fast.
- PT** Instruct the patient not to consume alcoholic beverages before the test.
- Draw the specimen before starting folate therapy.

## Abnormal findings

### ▲ Increased levels

Pernicious anemia  
Vegetarianism  
Recent blood transfusions

### ▼ Decreased levels

Folic acid deficiency anemia  
Hemolytic anemia  
Malnutrition  
Malabsorption syndrome  
(e.g., sprue, celiac disease)  
Malignancy  
Pregnancy  
Alcoholism  
Anorexia nervosa

notes

**fungal testing** (Antifungal antibodies; Beta-D-glucan,(1,3)- $\beta$ -D-glucan, Fungitell, Fungal culture, Fungal antigen assay, Fungal PCR testing)

**Type of test** Blood, microscopic examination

### Normal findings

No antibodies detected

$\beta$ -D-glucan:

Negative	< 60 pg/mL
Indeterminate	60-79 pg/mL
Positive	$\geq$ 80 pg/mL

Culture: no growth in 24 days

Gram stain: no fungus seen

### Test explanation and related physiology

Few fungal diseases can be diagnosed clinically; many are diagnosed by isolating and identifying the infecting fungus in the clinical laboratory. Fungal infections can be superficial, subcutaneous, or systemic (deep). The systemic fungal infections (mycoses) are the most important, for which serologic antibody testing is performed. In general, mycoses are caused by the inhalation of airborne fungal spores. In the United States, the most serious fungal infections are coccidioidomycosis, blastomycosis, histoplasmosis, and paracoccidioidomycosis. These infections start out as primary pulmonary infections. *Aspergillus*, *Candida*, and *Cryptococcus* systemic infections usually affect only those with compromised immunity.

Fungal antibody testing is not highly reliable. In general, this testing is used for screening for antibodies to dimorphic fungi (*Blastomyces*, *Coccidioides*, *Histoplasma* spp.) and the antigen of *Cryptococcus neoformans* during acute infection. Antibodies are present in only about 70% to 80% of infected patients. When positive, they merely indicate that the person has an active or has had a recent fungal infection. These antibodies can be identified in the blood or cerebrospinal fluid (CSF). Antibodies can be tested singularly or as a fungal panel.

(1,3)- $\beta$ -D-glucan is used to support the diagnosis of invasive fungal disease (IFD) in at-risk patients. Normally, serum contains low levels of (1,3)- $\beta$ -D-glucan, presumably from yeasts present in the alimentary and GI tract. (1,3)- $\beta$ -D-glucan is produced by most invasive fungal organisms. D-glucan becomes elevated well in advance of conventional clinical signs and symptoms of IFD.



As opportunistic infections, IFDs are common among patients with hematologic malignancies or AIDS. They account for a growing number of nosocomial infections, particularly among organ transplant recipients and other patients receiving immunosuppressive treatments. (1,3)- $\beta$ -D-glucan is produced by most invasive fungal organisms. *Blastomyces* and *Cryptococcus* produce very low levels of (1,3)- $\beta$ -D-glucan. Mucoromycetes do not produce (1,3)- $\beta$ -D-glucan. It is important to note that negative results do not exclude fungal etiology, especially in the early stages of infection.

Correlation of the patient clinical condition with culture results is necessary. Fungus can be cultured from blood, body fluids, CSF, fresh tissue, bronchopulmonary secretions, or swabs of the ear, nose, and throat or from urine. Accurate fungal culture is labor intensive and requires a highly experienced laboratory. Results are not available quickly.

### Interfering factors

- False-positive results can occur if a patient's intestinal tract is colonized with *Candida* spp.
- False-positive results occur in patients on hemodialysis using cellulose membranes.
- False-negative results occur in serum that is hemolyzed, icteric, lipemic, or turbid.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or serum separator
- Indicate on the laboratory slip the particular antibody or panel of antibodies that are to be tested.
- Because some patients with fungal infection may be immunocompromised, instruct them to check for signs of infection at the venipuncture site.

### Abnormal findings

#### ▲ Increased levels

Acute fungal infection

Previous systemic exposure to fungal disease

### notes

**galectin-3 (GAL-3)****Type of test** Blood**Normal findings** $\leq 22.1$  ng/mL**Test explanation and related physiology**

Heart failure progresses primarily by dilatation of the ventricular cardiac chamber through remodeling in fibrosis as a response to cardiac injury and/or overload. Galectin-3 (GAL-3) is a biomarker that appears to be actively involved in both the inflammatory and fibrotic pathways involved in remodeling. GAL-3 is a carbohydrate-binding lectin whose expression is associated with inflammatory cells, including macrophages, neutrophils, and mast cells. GAL-3 has been linked to cardiac remodeling in the setting of heart failure and a variety of other cardiac insults. Elevated levels are associated with increased risk of mortality.

**Interfering factors**

- Hemolysis increases GAL-3 levels.
- Heterophile antibodies (see p. 630) increase GAL-3 levels.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

**Abnormal findings****▲ Increased levels**

Congestive heart failure

**notes**

## gallbladder nuclear scanning (Hepatobiliary scintigraphy, Cholescintigraphy, DISIDA scanning, HIDA scanning)

**Type of test** Nuclear scan

### Normal findings

Gallbladder, common bile duct, and duodenum are visualized within 60 minutes after radionuclide injection. (This confirms patency of the cystic and common bile ducts.)

### Test explanation and related physiology

Through the use of iminodiacetic acid analogs (IDAs) labeled with technetium-99m ( $^{99m}\text{Tc}$ ), the biliary tract can be evaluated in a safe, accurate, and noninvasive manner.

Failure to visualize the gallbladder 60 to 120 minutes after injection of the radionuclide dye is virtually diagnostic of an obstruction of the cystic duct (acute cholecystitis). Delayed filling of the gallbladder is associated with chronic or acalculous cholecystitis. The identification of the radionuclide in the biliary tree, but not in the bowel, is diagnostic of common bile duct obstruction.

Gallbladder function can be numerically determined by calculating the capability of the gallbladder to eject its contents after the injection of a cholecystokinetic drug. It is believed that an ejection fraction lower than 35% indicates chronic cholecystitis or functional obstruction of the cystic duct. Ultrasound has largely replaced this test for the diagnosis of acute cholecystitis.

Occasionally, morphine sulfate is given intravenously during nuclear scanning. The morphine causes increased ampullary contraction. This reproduces the patient's symptoms of biliary colic and forces the bile containing the radionuclide into the gallbladder, shortening the expected time of visualization of the gallbladder.

### Contraindications

- Patients who are pregnant because of the risk of fetal damage

### Interfering factors

- If the patient has not eaten for more than 24 hours, the radionuclide may not fill the gallbladder. This would produce a false-positive result.
- ✦ The administration of opiates can prolong the time for gallbladder identification.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Assure the patient that he or she will not be exposed to large amounts of radioactivity.
- PT** Instruct the patient to fast for at least 2 to 4 hours before the test. This fasting is preferable but not mandatory.

### During

- Note the following procedural steps:
  1. After IV administration of a  $^{99m}\text{Tc}$ -labeled IDA analog (e.g., mebrofenin, disofenin), the right upper quadrant of the abdomen is scanned.
  2. Serial images are obtained over 1 hour.
  3. Subsequent images can be obtained at 15- to 30-minute intervals.
  4. If the gallbladder, common bile duct, or duodenum is not visualized within 60 minutes after injection, delayed images are obtained up to 4 hours later.
  5. When an *ejection fraction* is to be determined, the patient is given a fatty meal or cholecystokinin to evaluate emptying of the gallbladder. The gallbladder is continually scanned to measure the percentage of isotope ejected.
- Note that a radiologist performs this study in 1 to 4 hours in the nuclear medicine department.
- PT** Tell the patient that the only discomfort associated with this procedure is the intravenous (IV) injection of radionuclide.

### After

- Obtain a meal for the patient if indicated.

## Abnormal findings

Acute cholecystitis

Chronic cholecystitis

Acalculous cholecystitis

Common bile duct obstruction secondary to gallstones, tumor, or stricture

Cystic duct syndrome

## notes

## gallium scan

**Type of test** Nuclear scan

### Normal findings

Diffuse, low level of gallium uptake, especially in the liver and spleen

No increased gallium uptake within the body

### Test explanation and related physiology

A gallium scan of the total body may be performed 24, 48, and 72 hours after an IV injection of radioactive gallium. Most commonly, a single scan is performed 2 to 4 days after the gallium injection. Gallium is a radionuclide that is concentrated by areas of inflammation and infection, abscesses, and benign and malignant tumors. However, not all types of tumors will concentrate gallium. Lymphomas are particularly gallium avid. Other tumors that can be detected by a gallium scan include sarcomas, hepatomas, and carcinomas of the gastrointestinal (GI) tract, kidney, uterus, stomach, and testicle.

This test is useful in detecting metastatic tumor. However, to a large degree, positron emission tomography (PET) scans (see p. 719) have replaced gallium scans for the identification of malignancy. The gallium scan is useful in demonstrating a source of infection in patients with a fever of unknown origin. Gallium can be used to identify noninfectious inflammation within the body in patients who have an elevated sedimentation rate. PET scans are more commonly used to identify areas of acute infection.

Another method of scanning is called *single-photon emission computed tomography (SPECT) imaging*. This provides a more detailed image.

### Contraindications

- Patients who are pregnant unless the benefits outweigh the risks of fetal damage

### Interfering factors

- Recent barium studies will interfere with visualization of gallium within the abdomen.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.

## 434 gallium scan

- If ordered, administer a cathartic or enema to the patient to minimize increased gallium uptake in the bowel.

### During

- Note the following procedural steps:
    1. The unsedated patient is injected with gallium.
    2. A total-body scan may be performed 4 to 6 hours later by slowly passing a radionuclide detector over the body.
    3. Additional scans are usually taken 24, 48, and 72 hours later.
    4. During the scanning process, the patient is placed in the supine position and occasionally in the lateral position.
  - Note that a nuclear medicine technologist performs each scan in approximately 30 to 60 minutes. Repeated scanning is required. Repeated injections are not necessary.
- PT** Inform the patient that test results are interpreted by a physician trained in nuclear medicine and are usually available 72 hours after the injection.

### After

- PT** Assure the patient that only tracer doses of radioisotopes have been used and that no precautions against radioactive exposure to others are necessary.

### Abnormal findings

Tumor

Noninfectious inflammation

Infection

Abscess

### notes

## gamma-glutamyl transpeptidase (GGTP, $\gamma$ -GTP, Gamma-glutamyl transferase [GGT])

**Type of test** Blood

### Normal findings

Male and female age 45 years and older: 8-38 units/L or 8-38 international units/L (SI units)

Female younger than age 45 years: 5-27 units/L or 5-27 international units/L (SI units)

Elderly: slightly higher than adult level

Child: similar to adult level

Newborn: five times higher than adult level

### Test explanation and related physiology

This test is used to detect liver cell dysfunction, and it very accurately indicates even the slightest degree of cholestasis. This is the most sensitive liver enzyme in detecting biliary obstruction, cholangitis, or cholecystitis. As with leucine aminopeptidase and 5-nucleotidase (see pp. 561 and 649), the elevation of GGTP generally parallels that of alkaline phosphatase; however, GGTP is more sensitive. GGTP is not increased in bone diseases as is alkaline phosphatase. A normal GGTP level with an elevated alkaline phosphatase level implies skeletal disease. Elevated GGTP and alkaline phosphatase levels imply hepatobiliary disease. Another important clinical aspect of GGTP is that it can detect chronic alcohol ingestion. Therefore it is very useful in the screening and evaluation of patients with alcoholism. GGTP is elevated in approximately 75% of patients who chronically drink alcohol.

Why this enzyme is elevated after an acute myocardial infarction is not clear. It may represent the associated hepatic insult (if elevation occurs in the first 7 days) or the proliferation of capillary endothelial cells in the granulation tissue that replaces the infarcted myocardium. The elevation usually occurs 1 to 2 weeks after infarction.

### Interfering factors

- Values may be decreased in late pregnancy.
- ✚ Drugs that may cause *increased* GGTP levels include alcohol, phenobarbital, and phenytoin.
- ✚ Drugs that may cause *decreased* levels include clofibrate and oral contraceptives.

## 436 gamma-glutamyl transpeptidase

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red
- Patients with liver dysfunction often have prolonged clotting times.

### Abnormal findings

#### ▲ Increased levels

Hepatitis

Cirrhosis

Hepatic necrosis

Hepatic tumor or metastasis

Hepatotoxic drugs

Cholestasis

Jaundice

Myocardial infarction

Alcohol ingestion

Pancreatitis

Cancer of the pancreas

Epstein-Barr virus (infectious mononucleosis)

Cytomegalovirus infections

Reye syndrome

### notes



## gastric emptying scan

**Type of test** Nuclear scan

### Normal findings

Normal values are determined by type and quantity of radiolabeled ingested food.

Time	Lower normal limits (%)	Upper normal limits (%)
0 minutes		
30 minutes	70	
1 hour	30	90
2 hours		60
3 hours		30
4 hours		10

Values lower than normal represent abnormally fast gastric emptying. Values higher than upper limits represent delayed gastric emptying.

### Test explanation and related physiology

In this study, the patient ingests a solid or liquid “test meal” containing a radionuclide such as technetium (Tc). The stomach is then scanned until gastric emptying is complete. This study is used to assess the stomach’s ability to empty solids or liquids and to evaluate disorders that may cause a delay in gastric emptying, such as obstruction (caused by peptic ulcers or gastric malignancies) and gastroparesis. This scan is also useful in determining the rate of gastric emptying. This is helpful in the diagnosis of gastric retention secondary to gastroparesis or gastric obstruction. It is helpful in evaluating patients who have postcibal nausea, vomiting, bloating, early satiety, belching, or abdominal pain.

### Contraindications

- Patients who are pregnant or lactating, unless the benefits outweigh the risk of fetal or newborn injury

### Interfering factors

- ✦ Drugs that *decrease* gastric emptying time include anticholinergics, opiates, and sedative–hypnotics. These medications should be withheld for 2 days before testing.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Inform the patient that only a small dose of nuclear material is ingested. Reassure the patient that this is a safe dose.
- PT** Instruct the patient to keep on nothing by mouth (NPO) status after midnight on the day of the test.
- PT** Tell the patient that smoking is prohibited on the day of examination because tobacco can inhibit gastric emptying.

### During

- Note the following procedural steps:
  1. In the nuclear medicine department, the patient is asked to ingest a test meal. In the *solid-emptying* study, the patient eats scrambled egg whites containing Tc. In the *liquid-emptying* study, the patient drinks orange juice or water containing  $^{99m}\text{Tc}$  diethylenetriamine pentaacetic acid (DTPA) or indium-111 DTPA.
  2. After ingestion of the test meal, the patient is imaged by a gamma camera that records gastric images. Images are obtained for 2 minutes every 30 to 60 minutes until gastric emptying is complete. This may take several hours, although each particular timed scan takes only a few minutes.
  3. With the use of computer calculations of timed images, the rate of gastric emptying can be determined.

### After

- PT** Assure the patient that no radiation precautions need to be taken.

## Abnormal findings

Gastric obstruction caused by gastric ulcer or cancer

Nonfunctioning GI anastomosis

Gastroparesis

## gastrin

**Type of test** Blood

### Normal findings

Adult: 0-180 pg/mL or 0-180 ng/L (SI units)

Levels are higher in elderly patients.

Child: 0-125 pg/mL

### Test explanation and related physiology

Zollinger-Ellison (ZE) syndrome (gastrin-producing pancreatic tumor) and G-cell hyperplasia (overfunctioning of G cells in the distal stomach) are associated with high serum gastrin levels. Patients with these tumors have aggressive peptic ulcer disease. Unlike a patient with routine peptic ulcers, a patient with ZE syndrome or G-cell hyperplasia has a high incidence of complicated and recurrent peptic ulcers. It is important to identify this latter group of patients to institute more appropriate, aggressive medical and surgical therapy. The serum gastrin level is normal in patients with routine peptic ulcer and greatly elevated in patients with ZE syndrome or G-cell hyperplasia.

It is important to note that patients who are taking antacid peptic ulcer medicines or have had peptic ulcer surgery or have atrophic gastritis will have a high serum gastrin level. However, levels usually are not as high as in patients with ZE syndrome or G-cell hyperplasia.

Not all patients with ZE syndrome exhibit increased levels of serum gastrin. Some may have *top normal* gastrin levels, which makes these patients difficult to differentiate from patients with routine peptic ulcer disease. ZE syndrome or G-cell hyperplasia can be diagnosed in these top normal patients by *gastrin stimulation tests* with the use of calcium or secretin. Patients with these diseases have greatly increased serum gastrin levels associated with the infusion of these drugs.

### Interfering factors

- Peptic ulcer surgery creates a persistent alkaline environment, which is the strongest stimulant to gastrin.
- Ingestion of high-protein food can result in an increase in serum gastrin two to five times the normal level.
- ✘ Patients with diabetes who take insulin may have falsely *increased* levels.

## 440 gastrin

- ❖ Drugs that may *increase* serum gastrin levels include antacids and H<sub>2</sub>-blocking agents (e.g., esomeprazole, lansoprazole, omeprazole, pantoprazole, rabeprazole).
- ❖ Drugs that may *decrease* levels include anticholinergics and tricyclic antidepressants.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red
- PT Tell the patient to avoid alcohol for at least 24 hours.
- For the *calcium infusion test*, administer calcium gluconate intravenously. A preinfusion serum gastrin level is then compared with specimens taken every 30 minutes for 4 hours.
- For the *secretin test*, administer secretin intravenously. Preinjection and postinjection serum gastrin levels are taken at 15-minute intervals for 1 hour after injection.

### Abnormal findings

#### ▲ Increased levels

Zollinger-Ellison syndrome  
G-cell hyperplasia  
Pernicious anemia  
Atrophic gastritis  
Gastric carcinoma

Chronic renal failure  
Pyloric obstruction or  
gastric outlet obstruction  
Retained antrum after  
gastric surgery

### notes

## gastroesophageal reflux scan (GE reflux scan, Aspiration scan)

**Type of test** Nuclear scan

### Normal findings

No evidence of GE reflux

### Test explanation and related physiology

GE reflux scans are used to evaluate patients with symptoms of heartburn, regurgitation, vomiting, and dysphagia. They are also used to evaluate the medical or surgical treatment of patients with GE reflux. Finally, *aspiration scans* may be used to detect aspiration of gastric contents into the lungs.

### Contraindications

- Patients who cannot tolerate abdominal compression
- Patients who are pregnant or lactating unless the benefits outweigh the risks

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Assure the patient that no pain is associated with this test.
- PT** Instruct the patient not to eat anything after midnight.

#### During

- Note the following procedural steps:

##### *GE reflux scan*

1. The patient is placed in the supine position and asked to swallow a tracer cocktail (e.g., orange juice, diluted hydrochloric acid, and  $^{99m}\text{Tc}$  colloid).
2. Images are taken of the patient's esophageal area.
3. The patient is asked to assume other positions to determine whether GE reflux occurs and, if so, in what position.
4. A large abdominal binder that contains an air-inflatable cuff is placed on the patient's abdomen. This is insufflated to increase abdominal pressure.
5. Images are again taken over the esophageal area to determine whether any GE reflux occurs.

##### *Aspiration scans*

- Delayed images are made over the lung fields 24 hours after injection of technetium to detect esophago-tracheal aspiration.

## 442 gastroesophageal reflux scan

- Note that this procedure is performed in the nuclear medicine department in approximately 30 minutes.
- PT** Remind the patient that no discomfort is associated with this test.
- For infants being evaluated for achalasia, note that the tracer is added to the feeding or formula. Nuclear tracer films are then taken over the next hour, with 24-hour delayed films as needed.

### After

- PT** Assure the patient that he or she has ingested only a small dose of nuclear material. No radiation precautions need to be taken against the patient or his or her bodily secretions.

### Abnormal findings

GE reflux

Pulmonary aspiration

### notes

## gastrointestinal bleeding scan (Abdominal scintigraphy, GI scintigraphy)

**Type of test** Nuclear scan

### Normal findings

No collection of radionuclide in GI tract

### Test explanation and related physiology

The GI bleeding scan is a test used to localize the site of bleeding in patients who are having active GI hemorrhage. The scan also can be used in patients who have suspected intraabdominal hemorrhage from an unknown source.

Arteriography has limitations in its evaluation of GI bleeding. Arteriography can determine the site of bleeding only if the rate of bleeding exceeds 0.5 mL/min for detection. The GI bleeding scan has several advantages over arteriography. It can detect bleeding if the rate is greater than 0.05 mL/min. Also, with the use of  $^{99m}\text{Tc}$ -labeled red blood cells (RBCs), delayed films (as long as 24 hours) can be obtained, indicating the site of an intermittent or extremely slow intestinal bleed.

A GI scintigram is sensitive in locating the area of GI bleeding; however, it is not very specific in pinpointing the site or cause of bleeding. It is important to realize that this test can take 1 to 4 hours to obtain useful information. An unstable patient may need to go to surgery in minutes.

### Contraindications

- Patients who are pregnant or lactating, unless the benefits outweigh the risks
- Medically unstable patients

### Interfering factors

- Barium in the GI tract may mask a small source of bleeding.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Inform the patient that no pretest preparation is required.
- PT** Assure the patient that only a small amount of nuclear material will be administered.
- PT** Instruct the patient to notify the nuclear medicine technologist if he or she has a bowel movement during the test. Blood in the GI tract can act as a cathartic.

## 444 gastrointestinal bleeding scan

### During

- Note the following procedural steps:
    1. Ten millicuries of freshly prepared  $^{99m}\text{Tc}$ -labeled sulfur colloid is administered intravenously to the patient. If  $^{99m}\text{Tc}$ -labeled RBCs are to be used, 3 to 5 mL of the patient's own blood is combined with the  $^{99m}\text{Tc}$  and reinjected into the patient.
    2. Immediately after administration of the radionuclide, the patient is placed under a scintillation camera.
    3. Multiple images of the abdomen are obtained at short intervals (5-15 minutes).
    4. Detection of radionuclide in the abdomen indicates the site of bleeding. If no bleeding sites are noted in the first hour, the scan may be repeated at hourly intervals for as long as 24 hours.
  - Note that areas of the bowel hidden by the liver or spleen may not be adequately evaluated by this procedure. Also, the rectum cannot be easily evaluated because other pelvic structures (e.g., the bladder) obstruct the view.
  - Note that this test is usually performed in approximately 60 minutes by a technologist in nuclear medicine.
- PT** Tell the patient that the only discomfort associated with this study is the injection of the radioisotope.

### After

- PT** Assure the patient that only tracer doses of radioisotopes have been used and that radiation precautions are not necessary.

### Abnormal findings

Ulcer

Tumor

Angiodysplasia

Polyps

Diverticulosis

Inflammatory bowel disease

Aortoduodenal fistula

### notes



**genetic testing** (Breast cancer [*BRCA*] and ovarian cancer, Colon cancer, Cardiovascular disease, Tay–Sachs disease, Cystic fibrosis, Melanoma, Hemochromatosis, Thyroid cancer, Paternity [parentage analysis], Forensic genetic testing)

**Type of test** Blood; miscellaneous

### Normal findings

No genetic mutation

### Test explanation and related physiology

Genetic testing is used to identify a predisposition to disease, to establish the presence of a disease, to establish or refute paternity, or to provide forensic evidence used in criminal investigations. As research progresses and the Human Genome Project provides more information, precise and accurate methods of identification of normal and mutated genes are becoming more common.

Tests for defective genes known to be associated with certain diseases are now commonly used in screening people who have certain phenotypes and family histories compatible with a genetic mutation. Genetic testing is done in addition to a family history (pedigree).

### Breast cancer and ovarian cancer genetic testing

Inherited mutations in *BRCA* (BReast CAncer) genes indicate an increased susceptibility for development of breast cancer. The two genes in which mutations are most commonly seen are *BRCA1* and *BRCA2*. The *BRCA1* gene exists on chromosome 17, and *BRCA2* is on chromosome 13. These genes encode tumor suppressor proteins. More than half of women who inherit mutations will develop breast cancer by the age of 50 years compared with fewer than 2% of women without the genetic defect (Table 19). Also, see Table 20 for Interpretation of *BRCA* results.

The *BRCA* genes also confer an increased susceptibility for ovarian cancer. In the normal population, fewer than 2% of women develop ovarian cancer by age 70 years. Of women with mutations of the *BRCA1* gene, 44% develop ovarian cancer by that age. Ovarian cancer is less commonly associated with the *BRCA2* gene (20%). Furthermore, a woman with a *BRCA* mutation who has already had breast cancer has a 65% chance of developing a contralateral breast cancer in her lifetime (compared with < 15% of women without the genetic defect). A woman with breast cancer and a *BRCA* genetic defect has a 10 times greater risk of developing ovarian cancer as a second primary cancer compared with similar women without the mutated form of the gene.

**TABLE 19** Who should be tested for *BRCA* mutations?

Patient with breast cancer	Family history (with at least one characteristic)
Diagnosed at $\leq 50$ years of age	<ul style="list-style-type: none"> <li>No other family history</li> </ul>
Diagnosed with 2 primary breast cancers	
Triple negative breast cancer at $\leq 60$ years of age	
Breast cancer at any age	<ul style="list-style-type: none"> <li>One relative <math>\leq 50</math> years of age with breast cancer</li> </ul>
Diagnosed at any age	<ul style="list-style-type: none"> <li>One relative with ovarian cancer at any age</li> <li>Two relatives with breast cancer at any age</li> <li>Two relatives with pancreatic or prostate cancer</li> <li>Personal history of ovarian cancer</li> <li>Ashkenazi Jewish heritage</li> <li>First- or second-degree relative with <i>BRCA</i> mutation</li> <li>Close relative with male breast cancer</li> <li>Family history not needed.</li> </ul>
Male breast cancer at any age	
No personal history of breast cancer	Male breast cancer in a close relative

**TABLE 20** Interpretation of *BRCA* results

Result	Risk of breast cancer
Positive for a deleterious mutation	increased risk
Genetic variant; suspected deleterious	? increased risk
Genetic variant; favor polymorphism	Normal risk
Genetic variant; uncertain significance	? risk
No deleterious mutation	Normal risk

These mutations have an autosomal dominant inheritance pattern, indicating that women who inherit just one genetic defect can develop the phenotypic cancers. Men with *BRCA* genetic mutations (most commonly *BRCA2*) are at an increased risk for the development of breast, prostate, and colon cancer. In addition, they can pass the mutation to their daughters. Because *BRCA* is an autosomal dominant gene, 50% of the children are at risk.

### Colon cancer genetic testing

There are multiple forms of colon cancer strongly associated with family history. These genetic defects are inherited in an autosomal dominant fashion and are important for genome mismatch repair.

The most common form is familial adenomatous polyposis (FAP). The patient presents with over one hundred polyps in his or her colon—one or more of which may degenerate into cancer. FAP is caused by a genetic mutation in the 5q 21-22 (*APC*) gene on chromosome 5. These genes encode tumor suppression proteins.

Hereditary nonpolyposis colorectal cancer (HNPCC) syndrome is also known as Lynch syndrome. These patients are more difficult to recognize because they may have few polyps. HNPCC is associated most often with mutations (defective DNA mismatch repair) of *MLH1*, *MLH2*, and *MSH6*. HNPCC is associated with several other cancers (endometrial, gastric, and ovarian).

### Tay–Sachs disease genetic testing

Tay–Sachs disease is a GM2 gangliosidosis characterized by the onset of severe mental and developmental retardation in the first few months of life. Affected children become totally debilitated by 2 to 5 years of age and die by ages 5 to 8 years. Another form of the same disease is *late-onset Tay–Sachs* or chronic GM2, also known as gangliosidosis. The basic defect in affected patients is a mutation in the hexosaminidase A gene, which is on chromosome 15. This gene is responsible for the synthesis of hexosaminidase (HEX), an enzyme that normally breaks down a fatty substance called GM2 gangliosides. Ashkenazi (Eastern European) Jews and non-Jewish French Canadians, particularly those in the Cajun population in Louisiana, are affected most. This gene is inherited as an autosomal recessive gene. Carriers have one affected gene. Affected individuals have both defective genes. A *carrier couple* has a 25% chance of having a child affected with the disease.

At present, there is no treatment for the disease; it is important to identify carriers so that reproductive counseling can be provided. HEX A protein testing (see p. 500) has been extremely effective for identification of carriers and affected individuals. Both the test for the protein and that for the gene mutation are performed on a blood sample or on chorionic villus samples obtained during amniocentesis (see p. 49).

### **Cystic fibrosis genetic testing**

Cystic fibrosis (CF) is caused by a mutation in the cystic fibrosis transmembrane conductance regulator (*CFTR*) gene located on chromosome 7. This gene encodes the synthesis of a protein that serves as a channel through which chloride enters and leaves some types of epithelial cells. A mutation in this gene alters the cell's ability to regulate the chloride transport (see p. 233).

Currently, more than 1000 mutations can cause CF; several dozens of these account for 90% of the mutations in European Americans. However, the most common mutation, which accounts for 70% of the CF cases, is known as Delta F508.

CF is an autosomal recessive disease. A carrier has one mutated *CFTR* gene. The person affected by CF has mutations on both copies of *CFTR*. Genetic testing is now used to identify both carriers of CF and neonates with the disease, as well as to detect fetal disease during pregnancy. The sweat chloride test (see p. 857) is more easily performed and is a cheaper way to diagnose the disease in affected children. The use of genetic testing for CF is often limited to those with a family history of CF, partners of CF patients, and pregnant couples with a family history of CF. Genetic testing can be performed on blood samples or on samples taken during chorionic villus sampling (CVS, p. 241) or during amniocentesis (see p. 49).

### **Melanoma genetic testing**

Recent progress in the genetics of cutaneous melanoma has led to the identification of two melanoma susceptibility genes:

- The tumor suppressor gene *CDKN2A* encoding the p16 protein on chromosome 9 p21
- The *CDK4* gene on chromosome 12 q13

The p16 genetic mutation is by far the most common form of hereditary melanoma. Characteristics of familial melanoma include frequent multiple primary melanomas, early age of onset of first melanoma, and frequently the presence of atypical or

dysplastic nevi (moles). Family members with the following characteristics may consider testing for p16 genetic mutations:

- Multiple diagnoses of primary melanoma
- Two or more family members with melanoma
- Melanoma and pancreatic cancer
- Melanoma and a personal or family history of multiple atypical nevi
- Relatives of a patient with a p16 genetic mutation

The average age at diagnosis is 35 years for those with a mutation in p16 versus 57 years in the general population. P16 carriers also have an increased risk for pancreatic cancer.

### Hemochromatosis genetic testing

The diagnosis of hemochromatosis is traditionally made by using serum iron studies. When hereditary hemochromatosis (HH) is suspected, mutation analysis of the *hemochromatosis-associated HFE genes* (*C282Y* and *H63D*) is done. HH is an iron overload disorder that is considered to be the most common inherited disease in whites; it affects 1 in 500 individuals. Increased intestinal iron absorption and intracellular iron accumulation lead to progressive damage of the liver, heart, pancreas, joints, reproductive organs, and endocrine glands. Without therapy, men may develop symptoms between 40 and 60 years of age and women after menopause.

Patients with symptoms and early biochemical signs of iron overload consistent with HH should be tested. Relatives of individuals with HH should also be studied. HFE genotyping could improve disease outcomes of the disease. Serum iron markers are monitored at more frequent intervals if an HFE mutation is detected. Early initiation of phlebotomy therapy reduces the frequency or severity of hemochromatosis-related symptoms and organ damage.

### Thyroid cancer genetic testing

Genetic testing for RET germline mutations has shown 100% sensitivity and specificity for identifying those at risk for developing inherited medullary thyroid cancer (multiple endocrine neoplasia [MEN] 2A, MEN 2B, or familial medullary thyroid carcinoma [FMTC]).

Use of the genetic assay allows earlier and more definitive identification and clinical management of those with a risk for FMTC. FMTC is surgically curable if detected before it has spread to regional lymph nodes. Thus there is an emphasis on early detection and intervention in families that are affected by

## 450 genetic testing

MEN types 2A and 2B and FMTC, which account for one-fourth of medullary thyroid cancer cases.

### Cardiac genetic testing

Mutations in sarcomeric genes cause early onset cardiac channelopathies and cardiomyopathies. These are rare but potentially lethal heart conditions that include long QT syndrome (LQTS), catecholaminergic polymorphic ventricular tachycardia (CPVT), hypertrophic cardiomyopathy (HCM), arrhythmogenic right ventricular cardiomyopathy, and dilated cardiomyopathy (DCM). Patients with a sarcomeric gene mutation are nearly three times more likely to have an adverse cardiac outcome (cardiovascular death, nonfatal ischemic stroke, or progression to severe heart failure). Identifying patients with these genetic mutations can help diagnose a patient's disease, guide treatment options, and determine whether family members are at risk.

### Paternity genetic testing (parentage analysis)

DNA testing is the most accurate form of testing to prove or exclude paternity when the identity of the biological father of a child is in doubt. By comparing DNA variants in the mother and child, it is possible to determine variants that the child inherited from the biological mother. Thus any remaining DNA variation must have come from the biological father. If the DNA from the tested man is found to contain these paternal characteristics, then the probability of paternity can be determined. Testing is more than 99% accurate. Testing is so reliable that it is admissible in court. Testing can be done on a mouth swab or blood. Results are usually available in 1 to 3 weeks.

Many parents are given misinformation at the time of twin births as to whether the twins are identical or fraternal. DNA samples from siblings can be analyzed to indicate whether twins are identical or fraternal with an accuracy of greater than 99%.

Unfortunately, prenatal testing of the fetal components for paternity testing requires invasive testing such as CVS or amniocentesis. There are times, particularly in circumstances of rape, when early pregnancy paternity identification is desired. *Noninvasive prenatal paternity testing* can now be performed accurately by extracting and amplifying fetal chromosome alleles from maternal blood.

### Forensic genetic testing

Forensic DNA testing is used with increasing frequency in today's courtrooms because of its accuracy. In a courtroom, the reliability of the evidence can protect the individual and society as

a whole. Furthermore, DNA testing can be so conclusive that it often motivates plea bargaining and thereby reduces court time. It can quickly establish guilt or innocence beyond a reasonable doubt. Because DNA does not change and deteriorates very slowly even after death, testing can be performed on any body part, cadaver, or live person. Specimens considered adequate for DNA testing include blood, teeth, semen, saliva, bone, nails, skin scrapings, and hair. Forensic testing is also used for body identification.

### Contraindications

- Patients who are not emotionally able to deal with the results

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting is required.
- PT** Encourage all patients to receive genetic counseling.
- PT** Tell the patient the time it will take to have the results back.
- PT** Inform the patient of the high costs of genetic testing and that it may not be covered by all medical insurance plans.

#### During

- Obtain the specimen in a manner provided by the specialized testing laboratory.
  - **Blood:** A venous blood sample is collected in a lavender-top tube. Cord blood can be used for infants.
  - **Buccal swab:** A cotton swab is placed between the lower cheek and gums. It is twisted and then placed on a special paper or in a special container. Usually two to four swabs are requested.
  - **Amniotic fluid:** At least 20 mL of fluid is preferred.
  - **CVS:** 10 mg of cleaned villi are sent as prescribed by the testing laboratory.
  - **Product of conception:** 10 mg of placental tissue are preserved in a sterile medium.
  - **Other body parts:** As much tissue as is available is sent for testing.

#### After

- Apply pressure or a pressure dressing to the venipuncture site.

## 452 genetic testing

- Make sure that the patient has an appointment scheduled for obtaining the results. It is very upsetting for a patient and family to wait for the results.
- Arrangements should be made to ensure genetic and emotional counseling after results are obtained.

### **Abnormal findings**

Genetic carrier state

Affected state

### **Notes**



**genomic testing** (Oncotype DX, Genotyping, Genomics)

**Type of test** Microscopic examination

**Normal findings**

Low risk genomic score

**Test explanation and related physiology**

Genomic testing of cancers is a clinically validated, multigene assay that provides a quantitative assessment of the likelihood of the presence of a particular cancer or the risk for distant or local cancer recurrence. It can also assess the benefit of therapy in newly diagnosed cancer patients. In cancer treatment, the evaluation of the likelihood of distant recurrence is usually based on multiple pathologic factors, such as nodal status, tumor size, tumor grade, targeted receptors, and other pathologic findings. However, these factors may not be adequate to sufficiently quantify the recurrence risk to provide significant insight into the risks and benefits of therapy.

Genomic testing is designed to provide additional quantitative data regarding the patient's tumor biology to assist in clinical decision making regarding the use of anticancer therapies. By providing answers to key questions about the aggressiveness of and appropriate treatment for early-stage cancer, these practice-changing tests can help physicians and patients select the right treatment at the right time in each individual case. Genomic testing can help patients and their physicians optimize their cancer care and outcomes, enabling many patients to avoid unnecessary procedures and therapies, and saving the healthcare system billions of dollars in cost of unnecessary treatments. For example, a genomic test performed on the patient's tumor after part or all of the tumor has been excised can provide a risk score. If the patient's tumor genomics indicate a "low risk score," she or he has only a slight chance of the tumor acting aggressively and will only derive minimal or no benefit from aggressive anticancer therapy. Patients with tumor genomics indicating a "high risk" score may have a significant chance of the tumor acting aggressively and can experience considerable benefit from anticancer therapy.

At present, genomic testing is available for assessing cancers of the breast, bladder, lung, thyroid, prostate, and colon; lymphomas; and hematologic malignancies. Genomic testing is also being used to treat chronic benign diseases. This, along with molecular testing, offers the promise of very extensive

## 454 genomic testing

information regarding pathogenesis, prognosis, and treatment options for benign and malignant diseases.

Our use of genomic testing as it relates to cancer care refers to identifying particular panels of numerous genes of a tumor that statistically have been shown to provide more accurate prognostic and therapeutic data. Information gained from genomic testing is quite different from the information derived from other genetic testing, which identifies a gene causing a specific disease in a patient. For the sake of clarity, we have chosen to discuss genetic testing separately (see p. 445). Furthermore, genomic testing is different from other tumor markers determined by other molecular laboratory methods and may provide similar information regarding prognosis and treatment of cancer.

### Breast cancer genomic testing

The test is intended for use in all newly diagnosed patients with early-stage (stage I, II or IIIa), invasive breast cancer who have node-negative or node-positive (1-3), estrogen receptor-positive (ER+), HER2-negative disease. A score between 0 and 100 is provided. A low score means the cancer has a lower chance of returning and therefore the patient has a lower chance of benefiting from chemotherapy. A high score means the cancer has a higher chance of recurring, and the patient has a higher chance of benefiting from chemotherapy.

A second genomic test is available for noninvasive breast cancer to determine which women with this type of cancer are at high risk of having their cancer return, either as ductal carcinoma in situ (DCIS) or as an invasive carcinoma. The *Oncotype DX* breast cancer test for patients with DCIS provides an individualized prediction of the 10-year risk of local recurrence (either DCIS or invasive carcinoma). This diagnostic test helps guide treatment decision regarding extent of surgery, hormone therapy, and radiation therapy. *MammaPrint* and *Mammostrat* are other examples of genomic tests for breast cancer.

### Prostate cancer genomic testing

Genomic testing for prostate cancer was developed to help men with early-stage prostate cancer make the most informed treatment decision (specifically, aggressive therapy versus close surveillance). Men with genomic scores that are “low risk” have only a 3% chance that their disease will become life threatening and therefore may choose no treatment. Those with “high risk” scores would benefit from immediate surgery or radiation therapy. Genomic prostate tests build on traditional clinical

and pathologic factors to provide additional, clinically relevant insight into the underlying tumor biology. The result is a more precise and accurate assessment of risk, which helps more men avoid the lifelong complications associated with treatments they do not need, while directing aggressive therapy to men who require immediate treatment. Decipher, Prolaris, and ProMark are examples of genomic tests for prostate cancer.

### Colon cancer genomic testing

The genomic colon cancer test is used to predict recurrence risk in patients with stage II and stage III colon cancer to provide an individualized approach to treatment planning. The test can determine the likelihood that the cancer cells recur or spread after surgery. Based on a recurrence risk score, a better-informed decision can be made about whether chemotherapy is needed following surgery (stage II) or whether oxaliplatin should be added to the chemotherapy regimen after surgery (stage III). In this way, the test may help some patients avoid the complications of treatments they do not need, while directing others to the therapy that is most likely to benefit them.

Additionally, some genomic testing provides microsatellite instability (MSI) gene and mismatch repair (MMR) gene testing using immunohistochemistry to detect the protein expression of MSI and the 2 MMR genes (MLH1 and MSH2) (see colon cancer tumor analysis and genetic testing, p. 445). Used in conjunction with the *Colon Recurrence Score* test, MMR status offers deeper insight into the risk of recurrence. Approximately 15% of stage II colon cancer patients with MMR-deficient (MMR-D) tumors have a lower risk of recurrence and adjuvant chemotherapy may provide minimal, if any, clinical benefit.

Furthermore, genomic testing and molecular testing can be used to identify particular cancer markers (see cancer tumor marker, p. 194) to indicate associative inherited syndromes, such as FAP and HNPCC (Lynch syndrome).

### Lung cancer genomic testing

See lung cancer molecular testing, p. 588.

### Thyroid cancer genomic testing

Thyroid cancer genomics can identify a benign gene expression signature in thyroid nodules that are identified as “indeterminate” on needle aspiration biopsy of the thyroid. In these cases, diagnostic thyroidectomy would not be required. This testing can also guide treatment decision when a “suspicious for malignancy” result is obtained. A thyroid nodule is most commonly

## 456 genomic testing

evaluated for the possibility of malignancy by ultrasound (see p. 891), thyroid scanning (see p. 883), thyroid hormone testing (see p. 894), and thyroid aspiration (see p. 881). Genomic testing can be used on the aspirate (see thyroid fine needle biopsy, p. 881) of indeterminate nodules and classify them as either benign (< 5% risk for malignancy) or suspicious (> 50% risk for malignancy). With a benign score, observation or ultrasound follow-up could be recommended *in lieu* of thyroid surgery, avoiding unnecessary surgery. With a suspicious score, thyroid surgery is recommended.

Furthermore, if a thyroid aspirate is suspicious for cancer, genomic testing can also identify genes that can classify the tumor as medullary thyroid cancer (MTC). If positive, more aggressive lymph node surgery may be indicated. Thyroid genomics of a suspicious aspirate could also identify the V600E *BRAF* mutation indicating the presence of a papillary thyroid cancer (PTC) and encourage more extensive planned surgery. The presence of the *BRAF* gene mutation may be a prognostic marker of aggressive papillary thyroid cancer and is significantly associated with recurrence, lymph node metastases, extrathyroidal extension, and advanced stage in papillary thyroid cancer. Mutated *BRAF* is virtually absent in follicular, Hurthle cell, medullary carcinomas, and in benign thyroid tumors. The *Afirma gene expressive classifier (GEC)* is an example of genomic test for thyroid cancer.

### Lymphoma and hematologic cancer genomics

A major challenge in diagnosing hematologic and solid tumor cancers is the high degree of heterogeneity of the tumor cells. Mutations that have critical clinical implications regarding diagnosis and treatment may only be present in very low levels, making detection of these mutations difficult. Genomic testing can identify mutations that can provide insight into prognosis by providing a comprehensive view of the tumor's genomic profile. Particular testing panels have been developed for follicular and B cell lymphomas and for other non-Hodgkin lymphomas and leukemias. These tests are performed with histochemistry, karyotyping, fluorescence in situ hybridization, and other molecular testing.

### Bladder cancer genomics

This genomic assay is able to provide an assessment of bladder tumor aggressiveness for individual patients. It is able to separate less aggressive bladder cancer from the more aggressive muscle invasive bladder cancer. The test is also able to predict the response and survival outcome to neoadjuvant chemotherapy.

Another bladder cancer–related genomic assay is performed on urine. This test is designed to help rule out the presence of bladder (or any urothelial) cancer at an early stage. It is particularly appropriate for low-risk patients with hematuria. This is also a noninvasive and accurate method of performing bladder cancer surveillance because of its ability to detect early recurrence of bladder cancer. A “low risk” score would indicate no evidence of a new or recurrent urothelial cancer. The *Cx Bladder* is an example of a genomic test for bladder cancer.

## Procedure and patient care

### Before

- PT** Explain the significance of the prognostic data available for the patient’s tumor.
- PT** Explain the benefits of genomics in helping the physician and the patient make appropriate decisions regarding the use of adjuvant chemotherapy.
- Provide the patient with emotional support through the post-operative period.
  - Ensure that the patient’s insurance will cover this expensive testing.

### During

- The pathologist will send paraffin-embedded tumor tissue to the centralized laboratory.
- Results will be available in about 1 to 3 weeks.
- For bladder cancer identification, a urine specimen is provided.

### After

- PT** Provide education and support to patients as they evaluate their results.

## Abnormal findings

High risk for the presence of cancer

Aggressive versus nonaggressive cancer

High risk for cancer recurrence

High likelihood of benefit for aggressive anticancer treatment

## notes

## gliadin, endomysial, and tissue transglutaminase antibodies

**Type of test** Blood

### Normal findings

	Age	Normal
Gliadin IgA/IgG	0-2 years	< 20 EU
	3 years and older	< 25 EU
Endomysial IgA	All	Negative
Tissue transglutaminase IgA	All	< 20 EU

### Test explanation and related physiology

Gliadin and gluten are proteins found in wheat and wheat products. Patients with celiac disease cannot tolerate ingestion of these proteins or any products containing wheat. When an affected patient ingests wheat-containing foods, gluten and gliadin build up in the intestinal mucosa. These gliadin and gluten proteins (and their metabolites) cause direct mucosal damage.

Furthermore, IgA immunoglobulins (antigliadin, antientomysial, and antitissue transglutaminase [tTG-ab]) are made, appearing in the gut mucosa and in the serum of severely affected patients. The identification of these antibodies in the blood of patients with malabsorption is helpful in supporting the diagnosis of celiac sprue or dermatitis herpetiformis. However, a definitive diagnosis of celiac disease can be made only when a patient with malabsorption is found to have the pathologic intestinal lesions characteristic of celiac disease. Also, the patient's symptoms must be improved with a gluten-free diet. Both are needed for the diagnosis. Because of the high specificity of IgA endomysial antibodies for celiac disease, the test may obviate the need for multiple small bowel biopsies to verify the diagnosis.

In patients with known celiac disease, these antibodies can be used to monitor disease status and dietary compliance. Furthermore, these antibodies identify successful treatment because they will become negative in patients on a gluten-free diet.

### Interfering factors

- Other GI diseases (e.g., Crohn disease, colitis, and severe lactose intolerance) can cause elevated gliadin antibodies.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Obtain a list of foods that have been ingested in the last 48 hours.
- Assess how many malabsorption symptoms the patient has been experiencing in the last few weeks.

## Abnormal findings

Celiac disease

Celiac sprue

Dermatitis herpetiformis

notes

## glucagon

**Type of test** Blood

### Normal findings

50-100 pg/mL or 50-100 ng/L (SI units)

### Test explanation and related physiology

Glucagon is a hormone secreted by the alpha cells of the pancreatic islets of Langerhans. It is secreted in response to hypoglycemia and increases the blood glucose. As serum glucose levels rise in the blood, glucagon is inhibited by a negative feedback mechanism.

Elevated glucagon levels may indicate the diagnosis of a *glucagonoma* (i.e., an alpha islet cell neoplasm). Glucagon deficiency occurs with extensive pancreatic resection or with burned-out pancreatitis. Arginine is a potent stimulator of glucagon. If glucagon levels fail to rise even with arginine infusion, a diagnosis of glucagon deficiency as a result of pancreatic insufficiency is confirmed.

In an insulin-dependent patient with diabetes, glucagon stimulation caused by hypoglycemia does not occur. To differentiate the causes of glucagon insufficiency between pancreatic insufficiency and diabetes, *arginine stimulation* is performed. Patients with diabetes will have an exaggerated elevation of glucagon with arginine. In pancreatic insufficiency, glucagon is not stimulated with arginine. Furthermore, in patients with diabetes, hypoglycemia fails to stimulate glucagon release as would occur in a nondiabetic person.

Because glucagon is thought to be metabolized by the kidneys, renal failure is associated with high glucagon and, as a result, high glucose levels. When rejection of a transplanted kidney occurs, one of the first signs of rejection may be increased serum glucagon levels.

### Interfering factors

- Test results may be invalidated if a patient has undergone a radioactive scan within the previous 48 hours.
- Levels may be elevated after prolonged fasting or moderate to severe exercise.
- ✚ Drugs that may cause *increased* levels include some amino acids (e.g., arginine), danazol, gastrin, glucocorticoids, insulin, and nifedipine.
- ✚ Drugs that may cause *decreased* levels include atenolol, propranolol, and secretin.



## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: lavender

## Abnormal findings

### ▲ Increased levels

Familial hyperglucagonemia  
Glucagonoma  
Diabetes mellitus  
Chronic renal failure  
Severe stress including  
infection, burns, surgery,  
and acute hypoglycemia  
Acromegaly  
Hyperlipidemia  
Acute pancreatitis  
Pheochromocytoma

### ▼ Decreased levels

Idiopathic glucagon  
deficiency  
Cystic fibrosis  
Chronic pancreatitis  
Postpancreatectomy  
Cancer of the pancreas  
Diabetes mellitus

## notes

**glucose, blood** (Blood sugar, Fasting blood sugar [FBS])**Type of test** Blood**Normal findings**

Cord: 45-96 mg/dL or 2.5-5.3 mmol/L (SI units)

Premature infant: 20-60 mg/dL or 1.1-3.3 mmol/L

Neonate: 30-60 mg/dL or 1.7-3.3 mmol/L

Infant: 40-90 mg/dL or 2.2-5.0 mmol/L

Child &lt; 2 years: 60-100 mg/dL or 3.3-5.5 mmol/L

Child &gt; 2 years to adult:

Fasting: 70-110 mg/dL or &lt; 6.1 mmol/L (Fasting is defined as no caloric intake for at least 8 hours.)

Casual:  $\leq 200$  mg/dL (< 11.1 mmol/L) (Casual is defined as any time of day regardless of food intake.)

Adult: 74-106 mg/dL or 4.1-5.9 mmol/L

Elderly:

60-90 years: 82-115 mg/dL or 4.6-6.4 mmol/L

&gt; 90 years: 75-121 mg/dL or 4.2-6.7 mmol/L

**Possible critical values**

Adult male: &lt; 50 and &gt; 400 mg/dL

Adult female: &lt; 40 and &gt; 400 mg/dL

Infant: &lt; 40 mg/dL

Newborn: &lt; 30 and &gt; 300 mg/dL

**Test explanation and related physiology**

Through an elaborate feedback mechanism, glucose levels are controlled by insulin and glucagon. Serum glucose levels must be evaluated in relation to meal time. For example, a glucose level of 135 mg/dL may be abnormal if the patient is in the fasting state, but this level would be within normal limits if the patient had eaten a meal within the previous hour.

In general, true glucose elevations indicate diabetes mellitus; however, one must be aware of many other possible causes of hyperglycemia. Similarly, hypoglycemia has many causes. The most common cause is inadvertent insulin overdose in patients with brittle diabetes. If diabetes is suspected by elevated fasting blood levels, glycosylated hemoglobin (see p. 471) or glucose tolerance tests (see p. 467) can be performed.

Glycosylated hemoglobin (see p. 471) is now being performed more frequently to identify diabetes because this blood test represents blood sugar levels over the previous 120 days. Glucose determinations may be performed frequently in new

patients with insulin-dependent diabetes to monitor closely and adjust the insulin dosage. Finger stick blood glucose determinations are often performed before meals and at bedtime.

For minimally invasive monitoring, a small, sterile, disposable glucose-sensing device is inserted into the subcutaneous tissue (usually the arm). This sensor measures the change in glucose in the interstitial fluid. This information is recorded in a small beeper-sized monitor for 3 to 4 days. The monitor is taken to the doctor's office, where it is connected to a standard personal computer. Specialized software then downloads the stored information, and a more effective insulin regimen can be developed.

### Interfering factors

- Many forms of stress (e.g., general anesthesia, cerebrovascular accident, myocardial infarction, shock, strenuous exercise, burns) can cause increased serum glucose levels.
- Many pregnant women experience some degree of glucose intolerance. If significant, it is called gestational diabetes.
- ✚ Most IV fluids contain dextrose, which is quickly converted to glucose. Therefore most patients receiving IV fluids will have *increased* glucose levels.
- ✚ Drugs that may cause *increased* levels include antidepressants (tricyclics), antipsychotics, beta-adrenergic blocking agents, corticosteroids, cyclosporine, dextrose IV infusion, dextrothyroxine, diazoxide, diuretics, epinephrine, estrogens, glucagon, isoniazid, lithium, niacin, phenothiazines, phenytoin, salicylates (acute toxicity), statins, and triamterene.
- ✚ Drugs that may cause *decreased* levels include acetaminophen, alcohol, alpha-glucosidase inhibitors, anabolic steroids, biguanides, clofibrate, disopyramide, gemfibrozil, incretin mimetics, insulin, meglitinides, monoamine oxidase inhibitors, pentamidine, propranolol, sulfonylureas, and thiazolidinediones.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red or gray
- PT** For FBS, instruct the patient to fast for 8 hours. Water is permitted.
- PT** To prevent starvation, which may artificially raise the glucose levels, tell the patient not to fast much longer than 8 hours.
- Withhold insulin or oral hypoglycemics until after blood is obtained.

## 464 glucose, blood

- Glucose levels also can be evaluated by performing a finger stick and using either a *visually read test* or a *reflectance meter*. Using reflectance meters (e.g., glucometer, Accu-Chek bG, Stat-Tek) improves the accuracy of the blood glucose determination.

### Abnormal findings

#### ▲ Increased levels (hyperglycemia)

Diabetes mellitus  
Acute stress response  
Cushing syndrome  
Pheochromocytoma  
Chronic renal failure  
Glucagonoma  
Acute pancreatitis  
Diuretic therapy  
Corticosteroid therapy  
Acromegaly

#### ▼ Decreased levels (hypoglycemia)

Insulinoma  
Hypothyroidism  
Hypopituitarism  
Addison disease  
Extensive liver disease  
Insulin overdose  
Starvation

### notes

## glucose, postprandial (2-Hour postprandial glucose [2-Hour PPG], 1-Hour glucose screen)

**Type of test** Blood

### Normal findings

2-Hour PPG

0-50 years: < 140 mg/dL or < 7.8 mmol/L (SI units)

50-60 years: < 150 mg/dL

60 years and older: < 160 mg/dL

1-Hour glucose screen for gestational diabetes

< 140 mg/dL

### Test explanation and related physiology

For this study, a meal acts as a glucose challenge to the body's metabolism. Normally, insulin is secreted immediately after a meal in response to the elevated blood glucose level, causing the level to return to the premeal range within 2 hours. In patients with diabetes, the glucose level usually is still elevated 2 hours after the meal. The PPG is an easily performed screening test for diabetes mellitus. If the results are greater than 140 mg/dL and less than 200 mg/dL, a glucose tolerance test (GTT, p. 467) may be performed to confirm the diagnosis. If the 2-hour PPG is greater than 200 mg/dL, a diagnosis of diabetes mellitus is confirmed.

The *1-hour glucose screen* is used to detect gestational diabetes mellitus (GDM). Screening for GDM is performed with a 50- to 100-g oral glucose load followed by a glucose level determination 1 hour later. This is called the *O'Sullivan test*. Screening is done between weeks 24 and 28 of gestation. However, patients with risk factors such as a previous history of GDM may benefit from earlier screening. Patients whose serum glucose levels equal or exceed 140 mg/dL may be evaluated by a 3-hour glucose tolerance test. A 100-g oral glucose load is indicated for the diagnosis of gestational diabetes when results of the 50-g oral glucose load 1-hour screening test result are abnormal.

### Interfering factors

- Stress can increase glucose levels.
- If the patient is not able to eat the entire test meal or vomits some or all of the meal, levels will be falsely decreased.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.

## 466 glucose, postprandial

- Fasting: yes
- Blood tube commonly used: red or gray
- PT** For the *2-hour PPG*, instruct the patient to eat the entire meal (with at least 75 g of carbohydrates) and then not to eat anything else until the blood is drawn.
- For the *1-hour glucose screen for GDM*, give the fasting or non-fasting patient a 50-g oral glucose load.
- PT** Instruct the patient not to smoke during the testing. Smoking may increase glucose levels.
- PT** Inform the patient that he or she should rest during the 1- or 2-hour interval.

### Abnormal findings

#### ▲ Increased levels

Diabetes mellitus  
Gestational diabetes mellitus  
Malnutrition  
Hyperthyroidism  
Acute stress response  
Cushing syndrome  
Pheochromocytoma  
Chronic renal failure  
Glucagonoma  
Diuretic therapy  
Corticosteroid therapy  
Acromegaly  
Extensive liver disease

#### ▼ Decreased levels

Insulinoma  
Hypothyroidism  
Hypopituitarism  
Addison disease  
Insulin overdose  
Malabsorption or  
maldigestion

### notes

**glucose tolerance test (GTT, Oral glucose tolerance test [OGTT])**

**Type of test** Blood; urine

**Normal findings****Plasma test**

Fasting: < 110 mg/dL or < 6.1 mmol/L (SI units)

1 hour: < 180 mg/dL or < 10.0 mmol/L

2 hours: < 140 mg/dL or < 7.8 mmol/L

**Urine test**

Negative

**Test explanation and related physiology**


The GTT is used when diabetes is suspected (retinopathy, neuropathy, diabetic-type renal diseases). It is also suggested for the following:

- Patients with a family history of diabetes
- Patients who are markedly obese
- Patients with a history of recurrent infections
- Patients with delayed healing of wounds
- Women who have a history of delivering large babies, stillbirths, or neonatal births
- Patients who have transient glycosuria or hyperglycemia during pregnancy or after myocardial infarction, surgery, or stress

In the GTT, the patient's ability to tolerate a standard oral glucose load (usually 75 g of glucose) is evaluated by obtaining plasma and urine specimens for glucose level determinations before glucose administration and then at 1 hour and 2 hours afterward. Normally, there is a rapid insulin response to the ingestion of a large oral glucose load. This response peaks in 30 to 60 minutes and returns to normal in about 3 hours. Patients with an appropriate insulin response are able to tolerate the glucose load quite easily, with only a minimal and transient rise in plasma glucose levels within 1 to 2 hours after ingestion. In normal patients, glucose does not spill over into the urine.

Patients with diabetes will not be able to tolerate this load. As a result, their serum glucose levels will be greatly elevated from 1 to 5 hours (Figure 22). Also, glucose can be detected in their urine. It is important to note that intestinal absorption may vary among individuals. For this reason, some centers prefer the glucose load to be administered intravenously.

## 468 glucose tolerance test

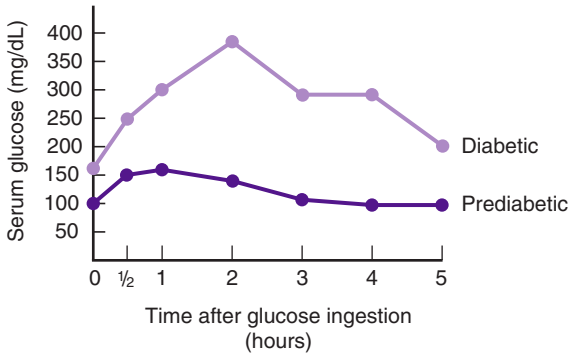

**FIG. 22** Glucose tolerance test curve for a diabetic and a prediabetic patient.

The American Diabetes Association recommends that pregnant women who have not previously had an abnormal GTT be tested at weeks 24 and 28 of gestation. A glucose level of more than 180 mg/100 mL 1 hour later is consistent with gestational diabetes.

Glucose intolerance also may exist in patients with oversecretion of hormones that have an ancillary effect on glucose, as in patients with Cushing syndrome, pheochromocytoma, acromegaly, aldosteronism, or hyperthyroidism. Patients with chronic renal failure, acute pancreatitis, myxedema, type IV lipoproteinemia, infection, or cirrhosis can also have an abnormal GTT. The GTT also is used to evaluate patients with reactive hypoglycemia. This may occur as late as 5 hours after the initial glucose load.

### Contraindications

- Patients with serious concurrent infections or endocrine disorders because glucose intolerance will be observed even though these patients may not have diabetes

### Potential complications

- Dizziness, tremors, anxiety, sweating, euphoria, or fainting during testing
- If these symptoms occur, a blood specimen is obtained. If the glucose level is too high, the test may need to be stopped and insulin administered.

### Interfering factors

- Smoking during the testing period stimulates glucose production because of the nicotine.



- Stress (e.g., from surgery, infection) can increase levels.
  - Exercise during the testing can affect glucose levels.
  - Fasting or reduced caloric intake before GTT can cause glucose intolerance.
- ✚ Drugs that may cause glucose intolerance include antihypertensives, antiinflammatory drugs, aspirin, beta-blockers, furosemide, nicotine, oral contraceptives, psychiatric drugs, steroids, and thiazide diuretics.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Educate the patient about the importance of having adequate food intake with adequate carbohydrates (150 g) for at least 3 days before the test.
- PT** Instruct the patient to fast for 12 hours before the test.
- PT** Instruct the patient to discontinue drugs (including tobacco) that could interfere with the test results.
- PT** Give the patient written instructions explaining the pretest dietary requirement.
- Obtain the patient's weight to determine the appropriate glucose loading dose (especially in children).

### During

- Obtain fasting blood and urine specimens.
- Administer the prescribed oral glucose solution, usually a 75- to 100-g carbohydrate load.
- Give pediatric patients a carbohydrate load based on their body weight.
- PT** Instruct the patient to ingest the entire glucose load.
- PT** Tell the patient that he or she cannot eat anything until the test is completed. However, encourage the patient to drink water.
- PT** Inform the patient that tobacco, coffee, and tea are not allowed because they cause physiologic stimulation.
- Collect a venous blood sample in a gray-top tube at 30 minutes and at hourly periods. Collect urine specimens at hourly periods.
- Mark on the tubes the time that the specimens are collected.
- Assess the patient for such reactions as dizziness, sweating, weakness, and giddiness. (These are usually transient.)
- For the IV GTT, administer the glucose load intravenously over 3 to 4 minutes.

### After

- Send all specimens promptly to the laboratory.
- Allow the patient to eat and drink normally.

## 470 glucose tolerance test

- Administer insulin or oral hypoglycemics if ordered.
- Apply pressure to the venipuncture site.

### Abnormal findings

Diabetes mellitus

Acute stress response

Cushing syndrome

Pheochromocytoma

Chronic renal failure

Glucagonoma

Acute pancreatitis

Diuretic therapy

Corticosteroid therapy

Acromegaly

Myxedema

Somogyi response to hypoglycemia

Postgastrectomy

notes

## glycosylated hemoglobin (GHb, GHB, Glycohemoglobin, Glycolated hemoglobin, Hemoglobin [Hb A<sub>1c</sub>], Diabetic control index)

**Type of test** Blood

### Normal findings

Nondiabetic adult or child: 4% to 5.9%

Good diabetic control: < 7%

Fair diabetic control: 8% to 9%

Poor diabetic control: > 9%

(Values vary with laboratory method used.)

### Test explanation and related physiology

This test is used to diagnose and monitor diabetes treatment. It measures the amount of hemoglobin A<sub>1c</sub> (HbA<sub>1c</sub>) in the blood and provides an accurate long-term index of the patient's average blood glucose level. In adults, about 98% of the hemoglobin in the RBCs is hemoglobin A. About 7% of hemoglobin A consists of a type of hemoglobin (HbA<sub>1</sub>) that can combine strongly with glucose in a process called *glycosylation*. When glycosylation occurs, it is not easily reversible.

HbA<sub>1</sub> is actually made up of three components: A<sub>1a</sub>, A<sub>1b</sub>, and A<sub>1c</sub>. HbA<sub>1c</sub> is the component that combines most strongly with glucose. Therefore HgA<sub>1c</sub> is the most accurate measurement because it contains the majority of glycosylated hemoglobin.

The amount of glycosylated hemoglobin (glycohemoglobin [GHb]) depends on the amount of glucose available in the bloodstream over a RBC's 120-day life span. Therefore determination of the GHb value reflects the average blood sugar level for the 100- to 120-day period before the test. The more glucose the RBC was exposed to, the greater the GHb percentage. One important advantage of this test is that the sample can be drawn at any time because it is not affected by short-term variations (e.g., food intake, exercise, stress, hypoglycemic agents).

As mentioned, the average life span of an RBC is 120 days, so the GHb may not reflect more recent changes in glucose levels. Because the turnover rate of proteins is much faster than hemoglobin, the measurement of serum *glycated proteins* (e.g., *glycated albumin* or *fructosamine*) provides more recent information about glucose levels. Glycated proteins reflect an average blood glucose level of the past 15 to 20 days. Although an initial single glycated protein result may not separate good glucose control from poor control, serial testing provides a much better indication of glucose control.

## 472 glycosylated hemoglobin

The GHb or glycated protein tests are particularly beneficial for the following:

- Evaluating patient compliance and success of treatment
- Comparing and contrasting the success of forms of diabetic therapy
- Determining the duration of hyperglycemia in patients with newly diagnosed diabetes
- Providing a sensitive estimate of glucose imbalance in patients with mild diabetes
- Individualizing diabetic control regimens
- Providing a feeling of reward for many patients when the test shows achievement of good diabetic control
- Evaluating the diabetic patient whose glucose levels change significantly day to day (brittle diabetic)
- Differentiating short-term hyperglycemia in patients who do not have diabetes (e.g., recent stress or myocardial infarction) from those who have diabetes (where the glucose has been persistently elevated)

A diagnosis of diabetes mellitus can be made on the basis of the results from two tests—fasting blood glucose (see p. 462) and GTT—performed on separate days that are close in time. By a relatively simple calculation, GHb can be accurately correlated with the daily *mean plasma glucose (MPG)* level, which is the average glucose level throughout the day. This has been very helpful for diabetics and health care professionals in determining and evaluating daily glucose goals. Each 1% change in GHb represents a change of approximately 35 mg/dL MPG. See [Table 21](#).

### Interfering factors

- Hemoglobinopathies can affect results because the quantity of hemoglobin A (and, as a result, HbA<sub>1c</sub>) varies considerably.

**TABLE 21** Correlation between glycosylated hemoglobin and mean plasma glucose

A <sub>1c</sub> (%)	Approximate MPG (mg/dL)	Interpretation
4	65	Nondiabetic range
5	100	Nondiabetic range
6	135	Nondiabetic range
7	170	ADA target
8	205	Action suggested

ADA, American Diabetes Association.

- Falsely elevated values occur when the RBC life span is lengthened.
- Abnormally low levels of proteins may falsely indicate normal glycosylated protein levels despite high glucose levels.
- Ascorbic acid may falsely indicate low levels of fructosamine.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: gray or lavender

### Abnormal findings

#### ▲ Increased levels

Newly diagnosed diabetic patient  
 Poorly controlled diabetic patient  
 Nondiabetic hyperglycemia (e.g., acute stress response, Cushing syndrome, pheochromocytoma, glucagonoma, corticosteroid therapy, or acromegaly)  
 Splenectomized patients  
 Pregnancy

#### ▼ Decreased levels

Hemolytic anemia  
 Chronic blood loss  
 Chronic renal failure

notes

## growth hormone (GH, Human growth hormone [HGH], Somatotropin hormone [SH])

**Type of test** Blood

### Normal findings

Men: < 5 ng/mL (mcg/L [SI units])

Women: < 10 ng/mL (mcg/L [SI units])

Children:

Newborn: 5-23 ng/mL (mcg/L [SI units])

1 week: 2-27 ng/mL (mcg/L [SI units])

1-12 months: 2-10 ng/mL (mcg/L [SI units])

1 year female: 0-10 ng/mL (mcg/L [SI units])

1 year male: 0-6 ng/mL (mcg/L [SI units])

### Test explanation and related physiology

This test is used to identify GH deficiency in adolescents who have short stature, delayed sexual maturity, or other growth deficiencies. It is also used to document the diagnosis of GH excess in patients with gigantism or acromegalia. GH is used to identify and follow patients with ectopic GH production by neoplasm. Finally, it is often used as a screening test for pituitary hypofunction or hyperfunction.

Because GH release is episodic, a random measurement of GH is unreliable to predict GH deficiency in adolescents. Measurement of free insulin-like growth factor (IGF) 1 and IGF BP 3 (see insulin-like growth factor, p. 541) is preferred in cases of short stature.

GH, or somatotropin, is secreted by the acidophilic cells in the anterior pituitary gland. It plays a central role in modulating growth from birth until the end of puberty. GH exerts its effects on many tissues through a group of peptides called *somatomedins*. The most commonly tested somatomedin is somatomedin C (also known as IGF-1), which is produced by the liver and has its major effect on cartilage.

If GH secretion is insufficient during childhood, limited growth and dwarfism may result. Also, a delay in sexual maturity may be a result in adolescents with reduced GH levels. Conversely, overproduction of GH during childhood results in gigantism, with the person sometimes reaching nearly 7 to 8 feet in height. An excess of GH during adulthood (after closure of long bone end plates) results in acromegaly, which is characterized by an increase in bone thickness and width but no increase in height.

Normal GH levels overlap significantly with deficient levels. Low GH levels may indicate deficiency or may be normal for certain individuals at certain times of the day. To negate time variables in GH testing, GH can be drawn 1 to 1½ hours after deep sleep has occurred. Levels increase during sleep. Also, strenuous exercise can be performed for 30 minutes in an effort to stimulate GH production.

To negate the common variations in GH secretion, screening for *IGF-1* or *somatomedin C* (see p. 541) provides a more accurate reflection of the mean plasma concentration of GH. These proteins are not affected by the time of day or food intake like GH is. A *GH stimulation test* (see p. 477) can be performed to evaluate the body's ability to produce GH. *Growth hormone suppression testing* is used to identify gigantism in children or acromegaly in the adult. If GH can be suppressed to < 2 ng/mL, neither of these conditions exists. The most commonly used suppression test is the oral glucose tolerance test (see p. 467). GH is normally suppressed when the glucose level increases. In acromegalic patients, only a slight decrease in GH occurs.

### Interfering factors

- Random measurements of GH are not adequate determinants of GH deficiency because hormone secretion is episodic.
- GH secretion is increased by stress, exercise, and low blood glucose levels.
- ✚ Drugs that may cause *increased* levels include amphetamines, arginine, dopamine, estrogens, glucagon, histamine, insulin, levodopa, methyl dopa, and nicotinic acid.
- ✚ Drugs that may cause *decreased* levels include corticosteroids and phenothiazines.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- The patient should not be emotionally or physically stressed because this can increase GH levels.
- It is preferred that the patient be fasting and well rested. Water is permitted.
- For *GH suppression testing*, the patient is kept NPO after midnight.

#### During

##### *Growth hormone test*

- Collect a venous blood sample in a red-top tube.

## 476 growth hormone

- Because approximately two-thirds of the total release of GH occurs during sleep, GH secretion also can be measured during hospitalization by obtaining blood samples when the patient is sleeping.

### *Growth hormone suppression test*

- Obtain peripheral venous access with normal saline solution.
- Obtain baseline GH and glucose levels as described previously.
- Administer the prescribed dose of glucose over 5 minutes.
- Obtain GH and glucose levels at 10, 60, and 120 minutes after glucose ingestion.

### **After**

- Indicate the patient's fasting status and the time the blood is collected on the laboratory slip. Include the patient's recent activity (e.g., sleeping, walking, eating).
- Send the blood to the laboratory immediately after collection.

### **Abnormal findings**

#### ▲ Increased levels

Gigantism  
Acromegaly  
Diabetes mellitus  
Anorexia nervosa  
Stress  
Major surgery  
Hypoglycemia  
Starvation  
Deep-sleep state  
Exercise

#### ▼ Decreased levels

Pituitary insufficiency  
Dwarfism  
Hyperglycemia  
Failure to thrive  
GH deficiency

### **notes**



## growth hormone stimulation test (GH provocation test, Insulin tolerance test [ITT], Arginine test)

**Type of test** Blood

### Normal findings

GH levels > 10 ng/mL or > 10 mcg/L (SI units)

IGF-1 > 80 ng/mL

### Test explanation and related physiology

Because GH (see p. 474) secretion is episodic, a random measurement of plasma GH is not adequate to make a diagnosis of GH deficiency. IGF-1 screening (see p. 541) followed by GH stimulation is indicated for children and adults suspected of GH deficiency. To diagnose GH deficiency, GH stimulation tests are sometimes needed. One of the most reliable GH stimulators is insulin-induced hypoglycemia, in which the blood glucose declines to less than 40 mg/dL. Other GH stimulants include vigorous exercise and drugs (e.g., arginine, clonidine, glucagon, levodopa). Glucagon is more widely used for GH stimulation because of safety concerns with insulin-induced hypoglycemia.

Usually a *double-stimulated test* is performed using an arginine infusion followed by insulin-induced hypoglycemia. A GH concentration of more than 10 mcg/L after stimulation effectively excludes the diagnosis of GH deficiency. Hypothyroidism should be excluded before GH stimulation testing.

### Contraindications

- Patients with epilepsy
- Patients with cerebrovascular disease
- Patients with myocardial infarction
- Patients with low basal plasma cortisol levels

### Potential complications

- Hypoglycemia so significant and severe as to cause ketosis, acidosis, and shock; with close observation, this is unlikely

### Procedure and patient care

#### Before

**PT** Explain the procedure very carefully to the patient and, if appropriate, to the parents.

**PT** Instruct the patient to remain NPO after midnight on the morning of the test. Water is permitted.

## 478 growth hormone stimulation test

### During

- Note the following procedural steps:
    1. A saline lock IV line is inserted for the administration of medications and for the withdrawal of frequent blood samples.
    2. Baseline blood levels are obtained for GH, glucose, and cortisol.
    3. Venous samples for GH are obtained at 60 and 90 minutes after injection of arginine, insulin, or glucagon.
    4. Blood glucose levels are monitored at 30-minute intervals with a glucometer. The blood sugar should drop to less than 40 mg/dL for effective measurement of GH reserve.
  - Monitor the patient for signs of hypoglycemia, postural hypotension, somnolence, diaphoresis, and nervousness. Ice chips are often given during the test for patient comfort.
  - This procedure is usually performed by a nurse with a physician in proximity.
  - This test takes approximately 2 hours to perform.
- PT** Tell the patient that the minor discomfort associated with this test results from the insertion of the IV line and the hypoglycemic response induced by the insulin injection.
- GH also can be stimulated by vigorous exercise. This entails running or stair climbing for 20 minutes. Blood samples of GH are obtained at 0, 20, and 40 minutes.

### After

- Observe the venipuncture site for bleeding.
  - Send the blood to the laboratory immediately after collection.
  - Give the patient cookies and punch or an IV glucose infusion.
- PT** Inform the patient and family that results may not be available for approximately 7 days. Some laboratories run GH tests only once per week.

### Abnormal findings

Growth hormone deficiency

Pituitary deficiency

### notes

## haptoglobin

**Type of test** Blood

### Normal findings

Adult: 50-220 mg/dL or 0.5-2.2 g/L (SI units)

Newborn: 0-10 mg/dL or 0-0.1 g/L (SI units)

### Possible critical values

< 40 mg/dL

### Test explanation and related physiology

The serum haptoglobin test is used to detect intravascular destruction (lysis) of red blood cells (RBCs), also called *hemolysis*. Haptoglobins, which are glycoproteins produced by the liver, are powerful, free hemoglobin (Hgb)-binding proteins. In hemolytic anemias associated with the hemolysis of RBCs, the released Hgb is quickly bound to haptoglobin, and the new complex is quickly catabolized. This results in a diminished amount of free haptoglobin in the serum; this decrease cannot be quickly compensated for by normal liver production. As a result, the patient demonstrates a transient, reduced level of haptoglobin in the serum. Megaloblastic anemias can reduce the haptoglobin level because of increased destruction of megaloblastic RBC precursors in the bone marrow.

Haptoglobins are also decreased in patients with primary liver disease not associated with hemolytic anemias. This occurs because the diseased liver is unable to produce these glycoproteins. Hematoma can reduce haptoglobin levels by the absorption of Hgb into the blood and binding with haptoglobin.

Elevated haptoglobin concentrations are found in many inflammatory diseases and can be used as a nonspecific acute-phase protein in much the same way as a sedimentation rate test is used (see p. 373). That is, levels of haptoglobin increase with severe infection, inflammation, tissue destruction, acute myocardial infarction, burns, and some cancers.

### Interfering factors

- A slight decrease in haptoglobin levels is noted in pregnancy.
- Ongoing infection can cause falsely elevated test results.
- ✎ Drugs that may cause *increased* haptoglobin levels include androgens and steroids.
- ✎ Drugs that may cause *decreased* levels include chlorpromazine, diphenhydramine, indomethacin, isoniazid, nitrofurantoin, oral contraceptives, quinidine, and streptomycin.

## 480 haptoglobin

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Collagen vascular disease  
Infection  
Tissue destruction  
Biliary obstruction  
Nephritis  
Pyelonephritis  
Ulcerative colitis  
Peptic ulcer  
Myocardial infarction  
Acute rheumatic disease  
Neoplasia

#### ▼ Decreased levels

Hemolytic anemia  
Transfusion reactions  
Prosthetic heart valves  
Systemic lupus erythematosus  
Primary liver disease not associated with hemolytic anemia  
Hemolytic disease of the newborn  
Hematoma  
Tissue hemorrhage  
Megaloblastic anemia  
Severe malnutrition

### notes

## Heinz body preparation

**Type of test** Blood

### Normal findings

No Heinz bodies detected

### Test explanation and related physiology

Heinz bodies are water-insoluble precipitates of oxidated–denatured Hgb that form within RBCs. They occur as a result of exposure to oxidative chemicals and drugs. Mutations of Hgb, thalassemias, and defects in the Hgb reductive defense system against oxidation lead to an enhanced tendency toward oxidative hemolysis. The diagnosis of these problems can be established by the detection of Heinz bodies in RBCs.

Heinz bodies are often associated with hemolytic anemias and the presence of spherocytosis.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender, pink, or green

### Abnormal findings

#### ▲ Increased levels

Unstable hemoglobinopathies (e.g., Hb Gun Hill)

RBC enzymopathies (e.g., G6PD)

Thalassemia

Heinz body hemolytic anemia

notes

***Helicobacter pylori* testing** (Anti-*Helicobacter pylori* antibody, *Campylobacter*-like organism [CLO] test, *H. pylori* stool antigen)

**Type of test** Blood, microscopic examination of antral or duodenal biopsy specimen, breath test, stool

**Normal findings**

**Blood test**

IgM	IgG
≤ 30 U/mL (negative)	< 0.75 (negative)
30.01-39.99 U/mL (equivocal)	0.75-0.99 (equivocal)
≥ 40 U/mL (positive)	≥ 1 (positive)

**Breath test**

No evidence of *H. pylori*

**Stool test**

No evidence of *H. pylori*

**Test explanation and related physiology**

*H. pylori*, a bacterium found in the mucus overlying the gastric mucosa and in the mucosa (cells that line the stomach), is a risk factor for gastric and duodenal ulcers, chronic gastritis, or even ulcerative esophagitis. There are several methods of detecting the presence of this organism. The organism can be cultured from a specimen of mucus obtained through a gastroscop (see p. 382). The organism can also be detected on a gastric mucosal biopsy (from the antrum and greater curvature of the corpus). This is very accurate. The gold standard for diagnosis of *H. pylori* disease is identifying the infected tissue by Gram, silver, Giemsa, or acridine orange stains.

It often takes several weeks before the results from cultures are available. It is preferable to start treatment before that time on a patient with symptomatic or active ulcer disease. For that reason, *rapid urease testing* for *H. pylori* was developed. *H. pylori* can break down large quantities of urea because of its ability to produce great amounts of an enzyme called urease. In the *CLO test*, a small piece of gastric mucosa (obtained through gastroscopy) is placed onto a specialized testing gel. If *H. pylori* organisms are present in the gastric mucosa, the urease (made by the *H. pylori*) will change the colors of the test material.

A *breath test* is also available for the detection of *H. pylori*. In the breath test, radioactive carbon ( $^{13}\text{C}$ ) is administered orally. The urea is absorbed through the gastric mucosa, where, if *H. pylori* is present, the urea will be converted to  $^{13}\text{CO}_2$  (where the carbon is radiolabeled). The  $^{13}\text{CO}_2$  is then taken up by the capillaries in the stomach wall and delivered to the lungs. There the  $^{13}\text{CO}_2$  is exhaled.

Although *H. pylori* does not survive in the stool, anti-*H. pylori* antibody can detect the presence of *H. pylori* antigen in a fresh stool specimen. Negative results indicate the absence of detectable antigen but do not eliminate the possibility of infection because of *H. pylori*.

*Serologic testing* is an inexpensive and noninvasive way of screening and diagnosing *H. pylori* infection. It is also used as a supportive diagnostic in which no preparation or abstinence from antacids is required. The IgG anti-*H. pylori* antibody is most commonly used. It becomes elevated 2 months after infection and stays elevated for more than 1 year after treatment. The IgA anti-*H. pylori* antibody, like IgG, becomes elevated 2 months after infection but decreases 3 to 4 weeks after treatment. The IgM anti-*H. pylori* antibody is the first to become elevated (about 3-4 weeks after infection) and is not detected 2 to 3 months after treatment. These antibody titers are fast becoming the gold standard for *H. pylori* detection. These antibodies can be detected with the use of a small amount of blood obtained by finger stick. Serologic testing is often used several months after treatment to document cure of *H. pylori* infection. Serologic testing is also used to corroborate the findings of other *H. pylori* testing methods.

### Contraindications

- Patients who are pregnant or are children
- The breath tests use radioactive carbon to which children should not be exposed.

### Interfering factors

- *H. pylori* can be transmitted by contaminated endoscopic equipment during endoscopic procedures.
- ✚ Rapid urease tests can be falsely negative if the patient uses antacid therapy within the week before testing.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

**PT** Tell the patient that no fasting is required for the blood test.

## 484 *Helicobacter pylori* testing

- If a biopsy or culture will be obtained by endoscopy, see discussion of esophagogastroduodenoscopy (EGD, p. 382).
- If culture is to be performed, be sure that the patient has not had any antibiotic, antacid, or bismuth treatment for 5 to 14 days before the endoscopy.

### During

- Collect a venous blood sample according to the protocol of the laboratory performing the test.
- A gastric or duodenal biopsy can be obtained by endoscopy.
- For the breath test, a dose of radioactive  $^{14}\text{C}$  or nonradioactive  $^{13}\text{C}$  urea is given by mouth.

### After

- Apply pressure to the venipuncture site.
- If endoscopy was used to obtain a culture, see procedure for EGD (see p. 382). The specimen should be transported to the laboratory within 30 minutes after collection.

## Abnormal findings

### ▲ Increased levels

- Acute and chronic gastritis
- Duodenal ulcer
- Gastric ulcer
- Gastric carcinoma

### notes



## hematocrit (Hct, Packed red blood cell volume, Packed cell volume [PCV])

**Type of test** Blood

### Normal findings

Male: 42%-52% or 0.42-0.52 volume fraction (SI units)

Female: 37%-47% or 0.37-0.47 volume fraction (SI units)

Pregnant female: > 33%

Elderly: values may be slightly decreased

Children (%)

Newborn: 44-64

2-8 weeks: 39-59

2-6 months: 35-50

6 months-1 year: 29-43

1-6 years: 30-40

6-18 years: 32-44

### Possible critical values

< 15% or > 60%


### Test explanation and related physiology

The Hct is a measure of the percentage of the total blood volume that is made up by the RBCs. The height of the RBC column is measured after centrifugation. It is compared with the height of the column of the total whole blood ([Figure 23](#)). The ratio of the height of the RBC column compared with the original total blood column is multiplied by 100%. This is the Hct value. It is routinely performed as part of a complete blood count. The Hct closely reflects the Hgb and RBC values. The Hct in percentage points usually is approximately three times the Hgb concentration in grams per deciliter when RBCs are of normal size and contain normal amounts of Hgb.

Abnormal values indicate the same pathologic states as abnormal RBC counts and Hgb concentrations (see p. 488). Decreased levels indicate anemia (reduced number of RBCs). Increased levels can indicate erythrocytosis. Like other RBC values, the Hct can be altered by many factors, such as hydration status and RBC morphology.

### Interfering factors

- Abnormalities in RBC size may alter Hct values.
- Extremely elevated white blood cell counts may affect values.
- Hemodilution and dehydration may affect the Hct level.

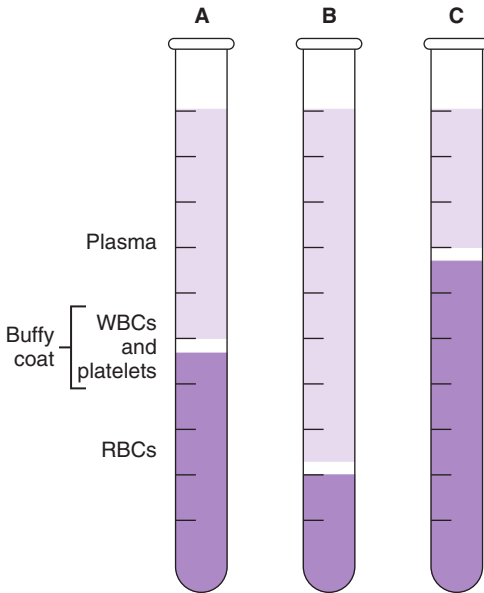

**FIG. 23** Tubes showing hematocrit levels of normal blood, blood with evidence of anemia, and blood with evidence of polycythemia. Note the buffy coat located between the packed red blood cells (RBCs) and the plasma. **A**, Normal percentage of RBCs. **B**, Anemia (low percentage of RBCs). **C**, Polycythemia (high percentage of RBCs). *WBC*, White blood cell.

- Pregnancy usually causes slightly decreased values because of hemodilution.
- Living in high altitudes causes increased values.
- Values may not be reliable immediately after hemorrhage.
- ✚ Drugs that may cause *decreased* levels include chloramphenicol and penicillin.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

**Abnormal findings****▲ Increased levels**

- Congenital heart disease
- Polycythemia vera
- Severe dehydration
- Erythrocytosis
- Eclampsia
- Burns
- Chronic obstructive pulmonary disease

**▼ Decreased levels**

- Anemia
- Hyperthyroidism
- Cirrhosis
- Hemolytic reaction
- Hemorrhage
- Dietary deficiency
- Bone marrow failure
- Normal pregnancy
- Rheumatoid arthritis
- Multiple myeloma
- Leukemia
- Hemoglobinopathy
- Prosthetic valves
- Renal disease
- Lymphoma
- Hodgkin disease

**notes**

**hemoglobin (Hb, Hgb)****Type of test** Blood**Normal findings**

Male: 14-18 g/dL or 8.7-11.2 mmol/L (SI units)

Female: 12-16 g/dL or 7.4-9.9 mmol/L (SI units)

Pregnant female: &gt; 11 g/dL

Elderly: values are slightly decreased

Children:

Newborn: 14-24 g/dL

0-2 weeks: 12-20 g/dL

2-6 months: 10-17 g/dL

6 months-1 year: 9.5-14 g/dL

1-6 years: 9.5-14 g/dL

6-18 years: 10-15.5 g/dL

**Possible critical values**

&lt; 5.0 g/dL or &gt; 20 g/dL

**Test explanation and related physiology**

The Hgb concentration is a measure of the total amount of Hgb in the peripheral blood, which reflects the number of RBCs in the blood. The test is normally performed as part of a complete blood count.

The Hct in percentage points usually is approximately three times the Hgb concentration in grams per deciliter when RBCs are of normal size and contain normal amounts of Hgb. Abnormal values indicate the same pathologic states as abnormal RBC counts and Hct concentrations (see p. 485). Decreased levels indicate anemia (reduced number of RBCs). Increased levels can indicate erythrocytosis. In addition, however, changes in plasma volume are more accurately reflected by the Hgb concentration. Slight decreases in the values of Hgb and the Hct during pregnancy reflect both the expanded blood volume caused by a chronic state of overhydration and an increased number of RBCs. Hemoglobinopathies, such as sickle cell disease and Hgb C disease, are also associated with reduced Hgb levels.

**Interfering factors**

- Slight Hgb decreases normally occur during pregnancy because of the expanded blood volume.
- Living in high-altitude areas causes high Hgb values.
- Heavy smokers have higher levels than nonsmokers.

- ▼ Drugs that may cause *increased* levels include gentamicin and methyldopa.
- ▼ Drugs that may cause *decreased* levels include antibiotics, anti-neoplastic drugs, aspirin, indomethacin, rifampin, and sulfonamides.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

### Abnormal findings

#### ▲ Increased levels

Congenital heart disease  
 Polycythemia vera  
 Hemoconcentration of the blood  
 Chronic obstructive pulmonary disease  
 Congestive heart failure  
 High altitudes  
 Severe burns  
 Dehydration

#### ▼ Decreased levels

Anemia  
 Hemorrhage  
 Hemolysis  
 Hemoglobinopathies  
 Nutritional deficiency  
 Lymphoma  
 Systemic lupus erythematosus  
 Sarcoidosis  
 Kidney disease  
 Chronic hemorrhage  
 Splenomegaly  
 Neoplasia

notes

## hemoglobin electrophoresis (Hgb electrophoresis)

**Type of test** Blood

### Normal findings

Adult/elderly: percentage of total Hgb

Hgb A<sub>1</sub>: 95%-98%

Hgb A<sub>2</sub>: 2%-3%

Hgb F: 0.8%-2%

Hgb S: 0%

Hgb C: 0%

Hgb E: 0%

Children: Hgb F

Newborn: 50%-80%

< 6 months: < 8%

> 6 months: 1%-2%

### Test explanation and related physiology

Hgb electrophoresis is a test that identifies and quantifies normal and abnormal forms of Hgb (hemoglobinopathies). Although many different forms of Hgb have been described, the more common types are A<sub>1</sub>, A<sub>2</sub>, F, S, E, and C. Each major Hgb type is electrically charged to varying degrees. Each electrophoretic band can be quantitated as a percentage of the total Hgb, indicating the severity of any recognized abnormality.

The form *Hgb A<sub>1</sub>* constitutes the major component of Hgb in the normal RBC. *Hgb A<sub>2</sub>* is only a minor component (2%-3%) of the normal Hgb total. *Hgb F* is the major Hgb component in a fetus but normally exists in only minimal quantities in a normal adult. Levels of Hgb F greater than 2% in patients older than age 3 years are considered abnormal. Hgb F is able to transport oxygen when only small amounts of oxygen are available (as in fetal life). In patients requiring compensation for prolonged chronic hypoxia (as in congenital cardiac abnormalities), Hgb F may be found in increased levels to assist in the transport of the available oxygen.

*Hgb S* and *Hgb C* are abnormal forms of Hgb that occur predominantly in African Americans. Hgb E occurs predominantly in Southeast Asians. The Hgb contents of some common disorders affecting Hgb, as determined by electrophoresis, are indicated in [Table 22](#) (see p. 491). Hgb E is produced less efficiently by RBC precursors; if there is an increased Hgb E content in the RBCs, those cells will have a low mean corpuscular volume (MCV, p. 772).

**TABLE 22** Hemoglobin (Hgb) contents of some common hemoglobinopathies

	Percentage range						
	Hgb A <sub>1</sub>	Hgb A <sub>2</sub>	Hgb F	Hgb S	Hgb H	Hgb C	Hgb E
Sickle cell disease	0	2-3	2	95-98	0	0	0
Sickle cell trait	50-65	2-3	2	35-45	0	0	0
Hgb C disease	0	2-3	2	0	0	90-100	0
Three gene deletion $\alpha$ -thalassemia (Hgb H disease)	65-90	0.3-1.5	0.6-4.5	0	0-30	0	0
$\beta$ -Thalassemia major	0	0-15	85-100	0	0	0	0
$\beta$ -Thalassemia trait	50-85	4-8	1-5	0	0	0	0
Hgb E disease	0	0	0	0	0	0	100

## 492 hemoglobin electrophoresis

Quantifying abnormal hemoglobins is helpful in determining the zygosity of a familial hemoglobinopathy. Furthermore, quantification of abnormal Hgb proteins provides a method of monitoring treatments designed to increase more effective Hgb variants and decrease abnormal variants.

### Interfering factors

- Blood transfusions within the previous 12 weeks may alter test results.
- Glycosylated Hgb can blur the peak of Hgb F and cause falsely low levels of Hgb F.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

### Abnormal findings

Sickle cell disease

Sickle cell trait

Hgb C disease

Hgb H disease

Thalassemia major

Thalassemia minor

Hgb E disease

notes



## hepatitis virus studies

**Type of test** Blood

### Normal findings

Negative

### Test explanation and related physiology

Hepatitis is an inflammation of the liver caused by viruses, alcohol ingestion, drugs, toxins, or overwhelming bacterial sepsis. The three common viruses now recognized to cause disease are hepatitis A, hepatitis B, and hepatitis C (also called non-A/non-B) viruses. Hepatitis D and E viruses are much less common in the United States. They are all associated with elevations of hepatocellular enzymes, such as aspartate aminotransferase (AST) and alanine aminotransferase (ALT).

*Hepatitis A virus (HAV)* was originally called *infectious hepatitis*. It has a short incubation period of 2 to 6 weeks and is highly contagious. During active infection, HAV is excreted in the stool and transmitted via oral–fecal contamination of food and drink. Most infections are not associated with symptoms severe enough to warrant medical evaluation. IgG and IgM antibodies to HAV are routinely used when HAV infection is suspected.

The first HAV antibody to appear is the IgM antibody (*HAV-Ab/IgM*) in approximately 3 to 4 weeks after exposure or just before hepatocellular enzyme elevations occur. These IgM levels usually return to normal in approximately 8 weeks. The next HAV antibody to rise is IgG (*HAV-Ab/IgG*), which appears approximately 2 weeks after the IgM begins to increase and slowly returns to normal levels. The IgG antibody can remain detectable for more than 10 years after the infection. If the IgM antibody is elevated in the absence of the IgG antibody, acute hepatitis is suspected. If, however, IgG is elevated in the absence of IgM elevation, a convalescent or chronic stage of HAV viral infection is indicated.

These antibodies may not be positive soon after infection occurs, which delays the investigation of infectious outbreaks. The HAV virus can be detected directly by measuring *HAV RNA* in the sera of patients suspected of acute infection.

*Hepatitis B virus (HBV)* is commonly known as *serum hepatitis*. It has a long incubation period of 5 weeks to 6 months. HBV is most frequently transmitted by blood transfusion; however, it also can be contracted via exposure to other body fluids. The incidence of hepatitis B is increased among blood

transfusion recipients, male homosexuals, dialysis patients, transplant patients, IV drug abusers, and patients with leukemia or lymphoma. Hospital personnel are also at increased risk of infection mostly because of needlestick contamination.

HBV, also called the *Dane particle*, is made up of an inner core surrounded by an outer capsule. The outer capsule contains the *hepatitis B surface antigen (HBsAg)*, formerly called *Australian antigen*. The inner core contains *HBV core antigen (HBcAg)*. The *hepatitis B e-antigen (HBeAg)* is also found in the core. Antibodies to these antigens are called HBsAb, HBcAb, and HBeAb. The tests used to detect these antigens and antibodies include the following (Table 23):

- *Hepatitis B surface antigen (HBsAg)*. This is the most frequently and easily performed test for hepatitis B, and it is the first test result to become abnormal. HBsAg rises before the onset of clinical symptoms, peaks during the first week of symptoms, and returns to normal by the time jaundice subsides. HBsAg generally indicates active infection by HBV. If the level of this antigen persists in the blood, the patient is considered to be a carrier or have chronic active hepatitis.
- *Hepatitis B surface antibody (HBsAb)*. This antibody appears approximately 4 weeks after the disappearance of the surface antigen and signifies the end of the acute infection phase. HBsAb also signifies immunity to subsequent infection. Concentrated forms of this agent constitute the hyperimmunoglobulin given to patients who have come in contact with HBV-infected patients. HBsAb is the antibody that denotes immunity after administration of hepatitis B vaccine.
- *Hepatitis B core antigen (HBcAg)*. No tests are currently available to detect this antigen.
- *Hepatitis B core antibody (HBcAb)*. This antibody appears approximately 1 month after infection with HBsAg and declines (although it remains elevated) over several years. HBcAb is also present in patients with chronic hepatitis. The HBcAb level is elevated during the time lag between the disappearance of HBsAg and the appearance of HBsAb. This interval is called the *core window*. During the core window, HBcAb is the only detectable marker of a recent hepatitis infection.
- *Hepatitis B e-antigen (HBeAg)*. This antigen generally is not used for diagnostic purposes but rather as an index of infectivity. The presence of HBeAg correlates with early

**TABLE 23** Hepatitis testing

Serologic findings	Appearance/disappearance	Application
HAV-Ab/IgM	4-6 weeks/3-4 months	Acute HAV infection
HAV-Ab/IgG	8-12 weeks/10 years	Previous HAV exposure or immunity
HAV RNA	First week of infection	Acute infection or carrier state
HBeAg	1-3 weeks/6-8 weeks	Acute HBV infection or chronic active hepatitis
HBeAb	4-6 weeks/4-6 years	Acute HBV infection ended, precore, or /core promoter mutant chronic infection
HBsAg	4-12 weeks/1-3 months	Acute or chronic HBV infection
HBsAb total	3-10 months/6-10 years	Previous HBV infection or immunity indicated
HBVc-Ab/IgM	2-12 weeks/3-6 months	Acute HBV infection
HBVc-Ab total	3-12 weeks/life	Previous HBV infection or convalescent stage
HBV DNA	First week of infection	Acute or chronic infection
HCV-Ab/IgG	3-4 months/2 years	Previous HCV infection
HCV RNA	First week of infection	Acute or chronic infection
HDV-Ag	1-3 days/3-5 days	Acute HDV infection or chronic infection
HDV-Ab/IgM	10 days/1-3 months	Acute HDV infection
HDV-Ab total	2-3 months/7-14 months	Chronic HDV infection

and active disease, as well as with high infectivity in acute HBV infection. The persistent presence of HBeAg in the blood predicts the development of chronic HBV infection.

- *Hepatitis B e-antibody (HBeAb)*. This antibody indicates that an acute phase of HBV infection is over, or almost over, and that infectivity is greatly reduced.

*Hepatitis B DNA* can be quantified and is a direct measurement of the *HBV viral load*. A one- or two-log decrease in viral load in a hepatitis B–infected patient means that antiviral therapy is working. A one- or two-log increase in a similar patient means that an antiviral has stopped working and that viral resistance may have developed. High levels of HBV DNA, ranging from 100,000 to more than 1 billion viral copies per milliliter, indicate a high rate of HBV replication. Low or undetectable levels, about 300 copies per milliliter or less, indicate an inactive infection. The World Health Organization established the international unit (IU) or copies per milliliter (mL), written as IU/mL or copies/mL, to measure HBV DNA.

*Hepatitis C (HCV)* (non-A/non-B [NANB] hepatitis) is transmitted in a manner similar to HBV. Most cases of hepatitis C are caused by blood transfusion. The incubation period is 2 to 12 weeks after exposure, and the clinical manifestations of the illness parallel those of HBV. However, unlike with HBV, HCV infection is chronic in more than 60% of infected persons. Although the disease course is variable, it is slowly progressive. Twenty percent of HCV patients develop cirrhosis and hepatocellular cancers associated with this chronic infection.

The screening test for detecting HCV infection is the detection of *anti-HCV antibodies* to HCV recombinant core antigen, NS3 gene, NS4 antigen, and NS5 antigen. The antibodies can be detected within 4 weeks of infection. With *HCV RNA testing*, the HCV virus can be directly detected and quantified. Like HBV DNA testing, HCV RNA viral load is usually expressed as units per milliliter or copies per milliliter. Although a higher viral load may not necessarily be a sign of more severe or more advanced disease, it does correlate with likelihood to respond to treatment. HCV RNA tests can also be used to monitor response to hepatitis C treatment.

*HCV genotypic testing*: The hepatitis C virus has seven different numbered DNA genotypes. Each of these genotypes has lettered subtypes. It may be important to find out the hepatitis C genotype because it could help determine both the type of treatment and the length of treatment; HCV genotype also helps to predict the likelihood of curing HCV. Worldwide and in

the United States, HCV genotype 1 is most common. It is possible to have more than one HCV genotype—this is more likely among injection drug users.

Hepatitis C home testing kits are now available using a small drop of blood obtained from a finger stick.

### Commonly used acute hepatitis testing

---

HAV-Ab IgM

HBsAg

HBcAb

HCV-Ab

*Hepatitis D virus (HDV)* is known to cause *delta hepatitis*. As stated earlier, HDV must enter the HBV to gain access to the liver and be infective. The patient must have HBV in the blood from a past or synchronously occurring infection. In the United States this is most commonly transmitted through tainted blood. The *HDV antigen* can be detected by immunoassay within a few days after infection. The IgM and total antibodies to HDV are also detected early in the disease. A persistent elevation of these antibodies indicates a chronic or carrier state.

*Hepatitis E virus (HEV)* was initially included in the non-A/non-B virus group but was isolated several years ago as an etiologic virus of short incubation. No antigen or antibody tests are currently available.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Usually, a hepatitis profile that includes several HBV antigens and antibodies is performed.

### Abnormal findings

#### ▲ Increased levels

Hepatitis A

Hepatitis B

Hepatitis C

Chronic carrier state, hepatitis B

Chronic hepatitis B

Hepatitis D

**herpes simplex** (Herpesvirus types 1 and 2, Herpes simplex virus types 1 and 2 [HSV 1, HSV 2], Herpes genitalis)

**Type of test** Blood; microscopic

**Normal findings**

No virus present

No HSV antibodies present

**Test explanation and related physiology**

HSV can be classified as either type 1 or type 2. *Type 1* is primarily responsible for oral lesions (blisters on the lips, or “cold sores”) or even corneal lesions. About half of the patients with HSV 1 develop recurrent infections. HSV 2 is a sexually transmitted viral infection of the urogenital tract.

Because most infants become infected if they pass through a birth canal containing HSV, determining its presence at delivery is necessary. Congenital infections may result in problems such as microcephaly, chorioretinitis, and mental retardation in the newborn. Disseminated neonatal herpes virus infections carry a high incidence of infant mortality. A vaginal delivery is possible if no virus is present, but birth by cesarean section is necessary if HSV is present.

Culture is still the standard criterion for HSV. Culture can be performed only during an outbreak. Serologic tests are more easily and conveniently available for detection of HSV 1 and HSV 2 antibodies. *Serologic tests for herpes simplex* are useful to supplement cultures or molecular detection for acute infection. The advantage of serology tests is that results can be available in a day. Serologic tests for IgG antibodies are available to help differentiate type 1 from type 2 infection. IgG antibodies indicate a previous exposure. IgM antibodies indicate an acute infection but do not differentiate well between types 1 and 2.

Fresh tissue is the definitive specimen for detection of HSV by molecular detection, particularly in infections involving oral, genital, central nervous system (CNS), ocular, and other sites.

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient.
- PT** Tell women to refrain from douching and tub bathing before the cervical culture is performed.
- Obtain the urethral specimen from men before voiding.

**During**

- Obtain cultures as follows:

**Urethral culture**

1. A culture is taken by inserting a sterile swab gently into the anterior urethra (see [Figure 39](#), p. 817) or genital skin lesion of the male patient.
2. Place the male patient in the supine position to prevent falling if vasovagal syncope occurs during introduction of the cotton swab or wire loop into the urethra.
3. The patient is observed for hypotension, bradycardia, pallor, sweating, nausea, and weakness.

**Cervical culture**

1. The female patient is placed in the lithotomy position, and a vaginal speculum is inserted.
2. Cervical mucus is removed with a cotton ball.
3. A sterile cotton-tipped swab is inserted into the endocervical canal and moved from side to side to obtain the culture. If a genital lesion is present, swabs from that area will be more sensitive in indicating infection.
4. For pregnant women with herpes genitalis, note that the cervix is cultured weekly for the herpes virus beginning 4 to 6 weeks before the due date. Vaginal delivery is possible if the following criteria are met:
  - a. The two most recent culture results are negative.
  - b. The woman is not experiencing any symptoms.
  - c. No lesions are visible on the vagina and vulva.
  - d. Throughout pregnancy, the woman has not had more than one positive culture, during which she was symptom-free.

**Blood for serology**

- Obtain a venous blood sample in a red-top tube.

**Molecular PCR tissue and other fluids**

- Obtain CSF (see p. 576) or other fluids by sterile technique as described elsewhere in this book.
- Obtain tissue by appropriate biopsy techniques.

**After**

**PT** Inform the patient how to obtain the test results.

**Abnormal findings**

Herpesvirus infection

notes

**hexosaminidase** (Hexosaminidase A, Hex A, Total hexosaminidase, Hexosaminidase A and B)**Type of test** Blood**Normal findings**

Hexosaminidase A: 7.5-9.8 units/L (SI units)

Total hexosaminidase: 9.9-15.9 units/L (SI units)

(Check with the laboratory because of the variety of testing methods.)

**Test explanation and related physiology**

Tay-Sachs disease (TSD) is a *lysosomal storage disease* (LSD), which, in infancy and early childhood, is characterized by loss of motor skills. Major categories of LSD include mucopolysaccharidoses, oligosaccharidoses, neuronal ceroid lipofuscinoses, and sphingolipidoses.

TSD, like other LSDs, is usually a result of a mutation in an autosomal recessive gene. Thus the affected person must have inherited a mutated gene from each parent in order to have TSD. Ashkenazi (Eastern European) Jews may particularly be carriers for these mutations. Eighty different mutations inhibit the function of this important gene. This gene encodes the synthesis of an enzyme called hexosaminidase. Without this enzyme, lysosomes of GM2 accumulate, particularly in the CNS.

Two clinically important isoenzymes of hexosaminidase have been detected in the serum: hexosaminidase A (hex A, made up of 1 alpha subunit and 1 beta subunit) and hexosaminidase B (hex B, made up of 2 beta subunits). Any genetic mutation that affects the alpha unit will cause a deficiency of hexosaminidase A, resulting in TSD. A mutation that affects the beta unit will cause a deficiency in hex A and B. Sandhoff disease, an uncommon variant of TSD, occurs with deficiency of both of these enzymes.

Biochemical testing to identify carriers and those affected by TSD is an important aspect of diagnosis. Hex A has been found to be abnormally low in carriers, whereas hex B is high. Therefore testing for total hexosaminidase is not useful. A carrier has a 25% chance of having a child with TSD if the other biological parent is also a carrier. In communities in which the Ashkenazi Jewish population is high, hex A screening has been very effective for identifying carriers. Furthermore, hex A is used to diagnose TSD in infants, young children, and adults.



Because accumulating analytes can be detected in urine, screening for LSDs typically begins with an analysis to detect disease-specific metabolite patterns. If positive, panel gene sequencing for LSDs (see p. 445) is useful to corroborate the identification of an affected person or a carrier.

### Interfering factors

- Hemolysis of the blood sample can cause inaccurate test results.
- Pregnancy can cause markedly increased values. For this reason, blood tests are not done during pregnancy.
- ▼ Oral contraceptives can falsely *increase* levels.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- PT Emphasize the importance of this test to Jewish couples of Eastern European ancestry who plan to have children. Explain that both must carry the defective gene to transmit TSD to their offspring.
- Professional genetic counseling should be provided to every person considering undergoing this test.
- Check with the laboratory regarding withholding contraceptives.
- Note that pregnant women can be evaluated by amniocentesis (see p. 49) or chorionic villus biopsy (see p. 241).
- Note that infants may have blood obtained by heel sticks. Neonates often have blood drawn through the umbilical cord.

### Abnormal findings

▼ Decreased  
hexosaminidase A  
Tay–Sachs disease

▼ Decreased  
hexosaminidase A and B  
Sandhoff disease

notes

## HIV drug resistance testing (HIV genotype)

**Type of test** Blood

### Normal findings

No resistant HIV

### Test explanation and related physiology

There are several factors that affect the success of HIV antiviral medications, including patient compliance, access to adequate care, optimal dosing, and drug pharmacology issues. Another significant factor that determines a patient's response to antiviral HIV drugs is the percentage of an HIV viral population that is resistant to the drugs that are administered. HIV resistance to therapy develops in 78% of patients. HIV genotyping is able to detect changes in the viral genome that are associated with drug resistance and is particularly able to predict HIV-1 resistance to protease and reverse-transcriptase inhibitor antiretroviral drugs.

HIV genotyping is particularly useful when failure of the most active antiviral therapy is suspected from decreasing CD4 counts (see p. 219). HIV genotyping can also be performed in conjunction with *HIV drug sensitivity testing*. HIV sensitivity testing estimates the ability of a cloned copy of the patient's virus to replicate in a cell culture in the presence of a particular antiviral drug. This same testing can help determine the amount of drug needed to inhibit viral replication.

### Interfering factors

- If the plasma HIV-1 RNA viral load is less than 1000 copies per mL of plasma, genotyping may be inaccurate.
- Minor HIV-1 populations that are less than approximately 20% of the total population may not be identified.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender or pink

**PT** Instruct the patient to observe the venipuncture site for infection. Patients with AIDS are immunocompromised and susceptible to infection.

### Abnormal findings

Drug resistance

notes

## HIV RNA quantification (HIV viral load)

**Type of test** Blood

### Normal findings

Undetected

### Test explanation and related physiology

Quantification of HIV RNA in the blood of patients infected with HIV can be used as an FDA-approved or supplementary test after serologic tests (see p. 506) are positive. Quantification is also helpful when confirmatory tests are indeterminate or cannot be accurately interpreted. Direct viral testing is helpful in differentiating newborn HIV infection from passive transmission of HIV antibodies from an HIV-infective mother. Finally, HIV RNA quantification testing determines HIV viral load. Determining viral load is used:

- to establish a baseline viral load before initiating anti-HIV-1 drug therapy
- to identify HIV-1 drug resistance while on anti-HIV therapy
- to identify noncompliance with anti-HIV-1 drug therapy
- to monitor HIV-1 disease progression
- to recommend the initiation of antiretroviral treatment (Table 24)
- to indicate the course of the disease because it is more accurate than any other test, including CD4 T-cell counts (see p. 219)
- as a determinant of patient survival (Table 25)

HIV viral load is most accurately determined by quantifying the amount of genetic material of the virus in the blood. In general, it is recommended to determine the baseline viral load by obtaining two measurements 2 to 4 weeks apart after HIV infection. Monitoring may continue with testing every 3 to 4 months in conjunction with CD4 counts. Both tests provide data used to determine when to start antiviral treatment. The viral load test can be repeated every 4 to 6 weeks after starting or changing antiviral therapy. Usually, antiviral treatment is continued until the HIV viral load is less than 500 copies/mL. It is important to recognize that a *nondetectable* result does not mean that no virus is left in the blood after treatment; it means that the viral load has fallen below the limit of detection by the test. A significant (> threefold) rise of viral load should warrant re-evaluation of therapy.

**TABLE 24** Recommendations for antiretroviral therapy based on viral load and CD4 count

CD 4 count ( $\times 10^5/L$ )	HIV RNA viral load, copies/mL		
	< 5000	5000-30,000	> 30,000
< 350	Recommend therapy	Recommend therapy	Recommend therapy
350-500	Consider therapy	Recommend therapy	Recommend therapy
> 500	Defer therapy	Consider therapy	Recommend therapy
Symptomatic		Recommend therapy	

**TABLE 25** Using the viral load to predict disease course

	HIV RNA viral load (copies/mL)				
	< 500	501-3000	3001-10,000	10,001-30,000	> 30,000
Developing AIDS (%)	5.4	16.6	31.7	55.2	80
Dying of AIDS (%)	0.9	6.3	18.1	34.9	69.5

## Interfering factors

- Incorrect handling and processing of the specimen can cause inconsistent results.
- Recent vaccinations may affect viral levels.
- Concurrent infections can cause inconsistent results.
- Variable compliance to therapy may alter test results.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender
- Specimens are often sent to a central laboratory.

**PT** Instruct the patient to observe the venipuncture site for infection. Patients with AIDS are immunocompromised and susceptible to infection.

**PT** Encourage the patient to discuss his or her concerns regarding the prognostic information from test results.

- Do not give test results over the phone. Increasing viral load results can have devastating consequences.
- Because test results vary according to the laboratory test method, it is important to use the same laboratory method for monitoring the course of the disease.
- Viral loads are usually repeated after starting or changing anti-viral therapy.

## Abnormal findings

HIV infection

notes

**HIV serologic and virologic testing** (AIDS serology, Acquired immunodeficiency serology, AIDS screen, Human immunodeficiency virus [HIV] antibody test, Western blot test, p24 direct antigen, HIV-RNA viral test)

**Type of test** Blood or fluid analysis (saliva)

### Normal findings

No evidence of HIV antigen or antibodies

### Test explanation and related physiology

There are two active types of human immunodeficiency viruses, types 1 and 2. HIV-1 is the most prevalent type within the United States and Western Europe, whereas HIV-2 is mostly limited to Western African nations. Serologic testing identifies antibodies developed as a result of HIV-1 or HIV-2 infections. Virologic tests identify RNA (or DNA) specific to HIV. Virologic tests can identify HIV infection in the first 11 days after infection. Serologic tests can identify HIV infection only after about 3 weeks. This 3-week time period is called the *seroconversion window*.

Serologic testing for HIV is divided into screening tests and confirmatory tests (Table 26).

In the past, serologic screening of patients suspected of having HIV-1 or HIV-2 infection usually began with an HIV antibody “screening test.” If positive, a confirmatory test was required to make the diagnosis of HIV infection. HIV serologic qualitative screening tests (for HIV-1 and HIV-2) were used to screen high- and low-risk individuals (Box 5) or for donor blood products. Because these rapid screening qualitative antibody immunoassays do not detect viral antigens, they could not detect infection in its earliest stage (before antibodies are formed). Because some persons who undergo HIV testing do not return to learn their test results, “point-of-service” rapid HIV antibody serologic screening testing in which results can be available in less than 1 hour is available. This is particularly helpful in urgent or emergent care points of service where HIV transmission could occur from blood or body fluid contamination. Furthermore, rapid antibody testing is helpful during labor in women whose HIV status is unknown. Point-of-care home kits are also available. Confirmatory tests for HIV-1 and HIV-2 antibodies include the Western blot assay and the indirect immunofluorescence assay (IFA). See screening recommendations (Table 27).

**TABLE 26** Serologic testing for HIV

Screening tests	Confirmatory or discriminatory tests
HIV-1 p24 antigen	WB HIV-1 antibody
HIV-1 antibody	WB HIV-2 antibody
HIV-2 antibody	Immunoblot–HIV-2 antibody
HIV-1/HIV-2 antibody	Immunofluorescence HIV-1 antibody (IFA)
Combined HIV-1/HIV-2 + HIV-1 p24 antigen	HIV RNA NAAT qualitative testing
Rapid HIV-1 antibody	
Rapid HIV-2 antibody	
Rapid HIV-1/HIV-2 antibody	

WB, Western blot; NAAT, nucleic acid amplification test.

#### BOX 5 Risk factors for HIV infection

- Sexually active male homosexuals
- Bisexual males
- Women with at-risk male partner
- Women with multiple male partners
- IV drug abusers
- Persons receiving blood products containing HIV
- Infants of HIV-positive mothers or mothers of unknown HIV status

**TABLE 27** Centers for Disease Control and Prevention HIV screening recommendations

Who	How often
All adults aged 18-64 years	Once in a lifetime
All adults with known risk factors	Yearly
All pregnant women	Once
Pregnant women at risk for HIV	Second test in third trimester
Newborns if mother is HIV positive or HIV status is unknown	Frequent repeated testing through first 6 months of life

## 508 HIV serologic and virologic testing

The *p24 direct serologic antigen assay* detects the viral protein p24 in the peripheral blood of HIV-infected individuals in which it exists either as a free (core) antigen or complexed to anti-p24 antibodies. The p24 antigen may be detectable as early as 16 days after infection. The p24 antigen test can be used to assess the antiviral activity of anti-HIV therapies. The p24 antigen test can also be used to differentiate active neonatal HIV infection from passive HIV antibody present from the mother's blood. It is also used to detect HIV infection before antibody seroconversion, to detect HIV in donor blood, and to monitor therapy.

The use of oral fluids for serologic HIV testing is as an alternative to serum testing. These HIV-1 antibody tests use *oral mucosal transudate (OMT)*, a serum-derived fluid that enters saliva from the gingival crevice and across oral mucosal surfaces. Another alternative to blood testing is *urine testing for HIV*. Only a spot urine collection is required. Testing urine for HIV antibodies is valuable, especially when venipuncture is inconvenient, difficult, or unacceptable. Insurance companies also commonly use it. It is important to note that all urine HIV tests are detecting antibodies and not the HIV particles. Urine does not contain the virus and is not a body fluid capable of infecting others. *HIV antigen/antibody (Ag/Ab) combination assays* are now available that can detect HIV infection on average 5 to 7 days earlier than assays that only detect antibodies.

### Interfering factors

- False-positive results can occur in patients who have autoimmune disease, lymphoproliferative disease, leukemia, lymphoma, syphilis, or alcoholism.
- False-positive results can occur in noninfected pregnant woman.
- HIV-2 infection can cause positive HIV-1/HIV-2 screening antibody test and an indeterminate WB HIV-1 confirmatory test.
- False-negative results can occur in the early incubation stage or end stage of AIDS.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Obtain an informed consent if required by the institution for all “opt in” testing.



- The blood is often sent to an outside laboratory for testing, although rapid antibody testing kits are becoming increasingly available in hospital laboratories and even in homes.
  - If the patient wishes to remain anonymous, use a number with the patient's name; be sure to record it accurately.
  - Note that if the serologic test is reactive (i.e., test result is positive twice consecutively), the Western blot test is performed on the same blood sample.
- PT** Instruct the patient to observe the venipuncture site for infection. Patients with AIDS are immunocompromised and susceptible to infection.
- Follow the institution's policy regarding test result reporting. Do not give results over the telephone. Remember that positive results may have devastating consequences, including loss of job, insurance, relationships, and housing.
- PT** Explain to the patient that a positive Western blot test merely implies HIV infection. It does not mean that the patient has clinical AIDS. Not all patients with positive antibodies will acquire the disease.
- PT** Encourage patients who test positive to inform their sexual contacts so that they can be tested. Most new cases of HIV infection are transmitted from patients who are unaware of their HIV status.
- PT** Inform the patient that subsequent sexual contact will put new partners at high risk for contracting AIDS.

## Abnormal findings

### ▲ Increased levels

AIDS

HIV infection

notes

## Holter monitoring

**Type of test** Electrodiagnostic


### Normal findings

Normal sinus rhythm

### Test explanation and related physiology

Holter monitoring is a continuous recording of the electrical activity of the heart. This can be performed for periods of up to 72 hours. With this technique, an electrocardiogram (ECG) is recorded continuously on magnetic tape during unrestricted activity, rest, and sleep (Figure 24). The Holter monitor is equipped with a clock that permits accurate time monitoring on the ECG tape. The patient is asked to carry a diary and to record daily activities as well as any cardiac symptoms that may develop during the period of monitoring.

Most units are equipped with an *event marker*. This is a button the patient can push when such symptoms as chest pain,

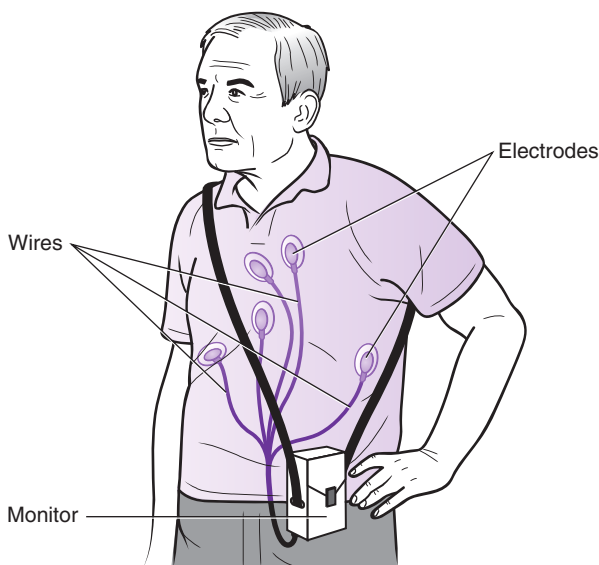

**FIG. 24** Electrical activity of the heart is recorded on a Holter monitor.

syncope, or palpitations are experienced. This type of monitor is referred to as an *event recorder*.

The Holter monitor is used primarily to identify suspected cardiac rhythm disturbances and to correlate these disturbances with symptoms (e.g., dizziness, syncope, palpitations, or chest pain). The monitor is also used to assess pacemaker function and the effectiveness of antiarrhythmic medications.

After completion of the determined time period, usually 24 to 72 hours, the Holter monitor is removed from the patient, and the record tape is played back at high speeds. The ECG tracing is usually interpreted by computer, which can detect any significant, abnormal waveform patterns that occurred during the testing.

*Implantable loop recorders (ILRs)* are used when long-term monitoring is required. These recorders are implanted subcutaneously via a small incision. ILRs can record ECG tracings continuously or only when purposefully activated by the patient. The ILR can be automatically activated by a predefined arrhythmia that will trigger device recording. If nothing irregular happens, then the information is subsequently erased. But if an arrhythmia does occur, the device locks it in and saves it to memory. ILRs can provide a diagnosis in many patients with unexplained syncope or presyncope.

### Contraindications

- Patients who are unable to cooperate with maintaining the lead placement from the monitor to the body
- Patients who are unable to maintain an accurate diary of significant activities or events

### Interfering factors

- Interruption in the electrode contact with the skin

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient about care of the Holter monitor.
- PT** Inform the patient about the necessity of ensuring good contact between the electrodes and the skin.
- PT** Teach the patient how to maintain an accurate diary. Stress the need to record significant symptoms.
- PT** Instruct the patient to note in the diary if any interruption in Holter monitoring occurs.
- PT** Assure the patient that the electrical flow is coming from the patient and that he or she will not experience any electrical stimulation from the machine.

## 512 Holter monitoring

- PT** Instruct the patient not to bathe during the period of cardiac monitoring.
- PT** Tell the patient to minimize the use of electrical devices (e.g., electric toothbrushes, shavers) that may cause artificial changes in the ECG tracing.

### During

- Prepare the sites for electrode placement with alcohol. (This is usually done in the cardiology department by a technologist.)
- Place the gel and electrodes at the appropriate sites. Usually, the chest and abdomen are the most appropriate locations for limb-lead electrode placement. The precordial leads also may be placed.
- PT** Encourage the patient to call if he or she has any difficulties.

### After

- Gently remove the tape and other paraphernalia securing the electrodes.
- Wipe the patient clean of electrode gel.
- PT** Inform the patient when the Holter monitoring interpretation will be available.

### Abnormal findings

Cardiac arrhythmia  
Ischemic changes

### notes

## homocysteine (Hcy)

**Type of test** Blood

### Normal findings

4-14  $\mu\text{mol/L}$

(Levels may increase with age.)

### Test explanation and related physiology

Homocysteine is an intermediate amino acid formed during the metabolism of methionine. Increasing evidence suggests that elevated blood levels of homocysteine may act as an independent risk factor for ischemic heart disease, cerebrovascular disease, peripheral arterial disease, and venous thrombosis. Homocysteine appears to promote the progression of atherosclerosis by causing endothelial damage, promoting low-density lipoprotein (LDL) deposition, and promoting vascular smooth muscle growth.

Dietary deficiency of vitamins B<sub>6</sub>, B<sub>12</sub>, or folate is the most common cause of elevated homocysteine. These vitamins are essential for the enzymatic metabolism of homocysteine to methionine (a protein). Homocysteine levels are elevated in patients with megaloblastic anemia. Some practitioners recommend homocysteine testing in patients with known poor nutritional status (alcoholics, drug abusers) and the elderly. Homocysteine is elevated in children with inborn errors of methionine metabolism.

Both fasting and postmethionine loading levels of homocysteine can be measured. In general, homocysteine levels less than 12 are considered optimal, levels from 12 to 15 are borderline, and levels greater than 15 are associated with high risk of vascular disease. When blood levels are elevated, urine homocysteine levels are also increased.

### Contraindications

- Patients whose creatinine levels exceed 1.5 mg/dL are not candidates for methionine loading. Elevated creatinine levels indicate malfunctioning kidneys that cannot effectively filter methionine.

### Interfering factors

- Patients with renal impairment have elevated levels of homocysteine because of poor excretion of the protein.
- Men usually have higher levels of homocysteine than women. Most likely, this is because of higher creatinine values and greater muscle mass. Values also increase with age.

## 514 homocysteine

- Smoking is associated with increased homocysteine levels.
- ✚ Drugs that may cause *increased* levels include azaribine, carbamazepine, methotrexate, nitrous oxide, and phenytoin.
- ✚ Oral contraceptives containing estrogen may alter the metabolism of homocysteine.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes (10-12 hours)
- Blood tube commonly used: purple, green, or blue
- For *methionine loading*, the patient ingests approximately 100 mg/kg of methionine after fasting for 10 to 12 hours. A blood sample is obtained. Repeat blood samples are collected at 2, 4, 8, 12, and 24 hours to compare levels of B vitamins and amino acids in the plasma.
- In the laboratory, the blood should be spun down within 30 minutes in order to avoid false elevation caused by release of homocysteine from RBCs.

### Abnormal findings

#### ▲ Increased levels

Cardiovascular disease  
Cerebrovascular disease  
Peripheral vascular disease  
Cystinuria  
Vitamin B<sub>6</sub> or B<sub>12</sub> deficiency  
Folate deficiency

### notes

## human chorionic gonadotropin (hCG, Pregnancy tests, hCG beta subunit)

**Type of test** Blood; urine

### Normal findings

- Negative: < 5 IU/L
- Indeterminate: 5-25 IU/L
- Positive: > 25 IU/L
- Males and nonpregnant females: < 2 IU/L

### Test explanation and related physiology

All pregnancy tests are based on the detection of human chorionic gonadotropin (hCG). hCG appears in the blood and urine of pregnant women within days after conception.

hCG is made up of alpha and beta subunits. The beta subunit is specific for hCG.

Very small levels of hCG can be detected and pregnancy can be determined 3 to 7 days after conception. There is no cross-over reactivity with other non-hCG glycoprotein hormones. The diagnostic cutoff for pregnancy is greater than 25 IU/L. Values between 5 and 25 IU/L are indeterminate for pregnancy. Results can be confirmed with a repeat test in 72 hours. Values in pregnancy should double every 3 days for the first 6 weeks. When an embryo is first large enough to be visible on transvaginal ultrasonography (see p. 685), the patient generally will have hCG concentrations between 1000 and 2000 IU/L. If the hCG value is high and gestational contents are not visible in the uterus, an ectopic pregnancy is suggested.

There are qualitative serum and urine hCG assays and quantitative serum hCG assays (Table 28). In the home setting, the urine is applied to a testing apparatus, and the color change is compared with a standard. If the color matches that standard,

**TABLE 28** Recommended uses for hCG testing

Test name	Recommended use
Qualitative beta hCG	Rapid pregnancy test
Quantitative hCG	More accurate pregnancy test
	Used to monitor high-risk pregnancy
Quantitative hCG (tumor marker)	Monitor patients with hCG-secreting tumors

## 516 human chorionic gonadotropin

pregnancy is present. Other test kits use the development of a line or plus symbol that may appear indicating pregnancy. These tests take only a few minutes to perform and to obtain results. They are best if performed a few days after a missed menses. However, results can be positive on the day of an expected menses.

Concentrations of hCG level off around week 20, significantly above prepregnancy levels. After delivery, miscarriage, or pregnancy termination, hCG falls until prepregnancy levels are reached.

Normally, hCG is not present in nonpregnant women. In a very small number of women (< 5%), hCG exists in minute levels. The presence of hCG does not necessarily indicate a normal pregnancy. Ectopic pregnancy, hydatidiform mole of the uterus, recent abortion, and choriocarcinoma can all produce hCG. However, hCG levels in ectopic pregnancy typically fail to double appropriately, and decreased levels eventually result relative to the values expected in normal intrauterine pregnancies of similar gestational age. hCG may be secreted by seminomatous and nonseminomatous testicular tumors, ovarian germ cell tumors, and benign or malignant nontesticular teratomas (see tumor markers, p. 194). Serial measurement of hCG after treatment is used to monitor therapeutic response in these tumors and will detect persistent or recurrent neoplastic disease.

### Interfering factors

- Tests performed too early in the pregnancy, before a significant hCG level exists, may cause false-negative results.
- Hematuria and proteinuria may cause false-positive results.
- Hemolysis of blood may interfere with test results.
- Urine pregnancy tests can vary according to the dilution of the urine. hCG levels may not be detectable in dilute urine but may be detectable in concentrated urine.
- ✔ Drugs that may cause false-negative urine results include diuretics (by causing diluted urine) and promethazine.
- ✔ Drugs that may cause false-positive results include anticonvulsants, antiparkinsonian drugs, hypnotics, and tranquilizers (especially promazine and its derivatives).

### Procedure and patient care

#### Before

- ✔ Explain the procedure to the patient.
- If a urine specimen will be collected, give the patient a urine container the evening before so that she can provide a first-voided morning specimen. This specimen generally contains the greatest concentration of hCG.



**During**

- Collect the first-voided urine specimen for urine testing.
- Collect a venous blood sample in a red-top tube for serum testing. Avoid hemolysis.

**After**

- Apply pressure to the venipuncture site.
- PT** Emphasize to the patient the importance of antepartal health care.

**Abnormal findings****▲ Increased levels**

Pregnancy  
Ectopic pregnancy  
Hydatidiform mole of the uterus  
Choriocarcinoma of the uterus, testes, or ovaries  
Tumor

**▼ Decreased levels**

Spontaneous abortion  
Fetal death  
Ectopic pregnancy

**notes**

**human lymphocyte antigen B27** (HLA-B27 antigen, Human leukocyte A antigen, White blood cell antigens, Histocompatibility leukocyte A antigen)

**Type of test** Blood

### Normal findings

Negative

### Test explanation and related physiology

The HLA antigens exist on the surface of white blood cells and on the surface of all nucleated cells in other tissues. The presence or absence of these antigens is determined genetically. Each gene controls the presence or absence of HLA A, B, C, or D.

The HLA system is used to assist in the diagnosis of certain diseases. For example, HLA-B27 is present in 80% of patients with Reiter syndrome. When a patient has recurrent and multiple arthritic complaints, the presence of HLA-B27 supports the diagnosis of Reiter syndrome. Other HLA–disease associations are mentioned in the abnormal findings section.

The HLA system of antigens has been used to indicate tissue compatibility in transplantation. Because HLA antigens are genetically determined, they can also be used to resolve *paternity investigations*.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: verify with laboratory

### Abnormal findings

#### ▲ Increased levels (HLA-B27 antigens present)

Ankylosing spondylitis  
 Reiter syndrome  
*Yersinia enterocolitica* arthritis  
 Anterior uveitis  
 Graves disease  
 Celiac disease or gluten enteropathy  
 Chronic active hepatitis  
 Multiple sclerosis  
 Myasthenia gravis  
 Dermatitis herpetiformis  
 Psoriasis

▲ **Increased levels (HLA-B27 antigens present)**

Juvenile diabetes

Hemochromatosis

Rheumatoid arthritis

notes

**human papillomavirus (HPV test, HPV DNA testing)****Type of test** Fluid analysis**Normal findings**

No HPV present

**Test explanation and related physiology**

An HPV test is performed to identify genital HPV infection in a woman with an abnormal PAP smear (see p. 667). HPV DNA incorporates itself into the cervical cell genome, promoting its effects through activation of oncogenes and suppression of host-cell immune response.

Genital HPV strains are divided into two groups (low, or nononcogenic risk, and high, or oncogenic risk), based on their oncogenic potential and ability to induce viral-associated tumors. Low-risk strains (HPV 6, 11, 42, 43, and 44) are associated with condylomata (genital warts) and low-grade cervical changes, such as mild dysplasia. Lesions caused by low-risk HPV infection have a high likelihood of regression and little potential for progression and are considered of no or low oncogenic risk. High-risk strains (HPV 16, 18, 31, 33, 35, 39, 45, 51, 52, 56, 58, 59, and 68) are associated with intraepithelial neoplasia and are more likely to progress to severe lesions and cervical cancer.

A clear causal relationship has been established between HPV infection and cervical cancer. Of the high-risk HPV strains, HPV 16 and 18 are the most carcinogenic and most prevalent. High-grade cervical intraepithelial lesions are most commonly associated with HPV 16 and 18, yet these strains are also frequently found to be the etiologic factor in minor lesions and mild dysplasia. Women who have normal Pap tests results and no HPV infection are at a very low risk (0.2%) for developing cervical cancer.

Numerous sources indicate that more than 60% of women with an abnormal Pap smear result will test positive for high-risk HPV. If the HPV test result is positive, the woman should undergo colposcopy to look for a more serious cervical lesion, such as cancer. It is well known that HPV infection in younger women is more prevalent and often spontaneously regresses, particularly in those younger than the age of 30 years. In contrast, persistent high-risk infection peaks in women older than 30 years. As a result, recent screening guidelines recommend that HPV testing be reserved for clinical use in the evaluation of women older than the age of 30 years and perhaps for younger women with high-grade squamous intraepithelial lesions. A combination

of cervical cytology and HPV DNA screening is also appropriate screening for women aged 30 years and older. Whether HPV testing can replace conventional Pap cytologic testing for cervical cancer screening awaits further study.

Testing for high-risk (oncogenic) HPV is summarized in [Table 29](#).

### Interfering factors

- Cervical specimens with low cellularity may diminish sensitivity of the test.
- High concentrations of antifungal cream or contraceptive jelly may diminish sensitivity of the test.

### Procedure and patient care

#### Before

**PT** Explain the procedure for Pap smear.

**PT** Instruct the patient not to douche or bathe in a tub during the 24 hours before the test.

**PT** Instruct the patient to empty her bladder.

**PT** Instruct the patient to reschedule testing if she is menstruating.

**PT** Tell the patient that no fasting or sedation is required.

#### During

- Note the following procedural steps:
  1. The patient is placed in the lithotomy position.
  2. With the use of either a cytology brush or a wooden spatula, a cervical mucous specimen is obtained by placing the instrument into the cervical os and rotating 3 to 5 times in clockwise and counterclockwise directions.

**TABLE 29** American Cancer Society recommendations for cervical screening

Population by age	Recommended screening
< 21 years	No screening
21-29 years	PAP or thin prep alone every 3 years
30-65 years	HPV and PAP/thin prep cotesting every 5 years
> 65 years	No screening after adequate negative prior screening
After hysterectomy	No screening
HPV vaccinated	Follow above recommendations (? delay screening for 3-5 years)

## 522 human papillomavirus

3. After specimen collection, rotate the broom-like device or spatula and Cytobrush several times in the collection vial to remove the specimen. Firmly cap the vial and discard the collection devices.
  4. Affix a patient identification label to the vial.
  5. Seal the vial and place in a plastic specimen bag along with a properly filled-out cytology requisition form and send them to the laboratory.
- Note that a smear is obtained by a nurse or a physician in approximately 10 minutes.
- PT** Tell the patient that no discomfort, except for insertion of the speculum, is associated with this procedure.

### After

- PT** Inform the patient that usually she will not be notified unless further evaluation is necessary.
- PT** Instruct the patient that HPV is a sexually transmitted disease. Proper precautions should be taken to prevent infecting sexual partners.

### Abnormal findings

HPV infection

notes

## human placental lactogen (hPL, Human chorionic somatomammotropin [HCS])

**Type of test** Blood

### Normal findings

Weeks of pregnancy	hPL concentration (mg/L = mcg/mL)
≤ 20	0.05-1
≤ 22	1.5-3
≤ 26	2.5-5
≤ 30	4-6.5
≤ 34	5-8
≤ 38	5.5-9.5
≤ 42	5-7

### Test explanation and related physiology

The human placenta produces hPL, which maintains pregnancy. Serum levels of hPL rise very early in normal pregnancy and continue to increase until a plateau is reached at about the 35th week after conception. Assays for maternal serum levels of hPL are useful in monitoring placental function. Measurements of hPL also are used in pregnancies complicated by hypertension, proteinuria, edema, postmaturity, placental insufficiency, or possible miscarriage.

A decreased serum concentration of hPL is pathognomonic for a malfunction of the placenta, which may cause intrauterine growth retardation, intrauterine death of the fetus, or imminent miscarriage. Pregnant women with hypertonia also show low serum concentrations of hPL. Because of the short biological half-life of hPL in serum, the determination of hPL gives a very accurate picture of the present situation.

Increased serum concentrations of hPL are found in women with diabetes mellitus and, because of the higher placental mass, in multiple pregnancies. In contrast to estriol, the hPL concentration depends only on the placental mass and not on fetal function. The simultaneous determination of hPL and estriol can be helpful in the differential evaluation of placental function.

### Interfering factors

- Prior nuclear medicine scans, because they can interfere with interpretation of test results

## 524 human placental lactogen

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: no
  - Blood tube commonly used: red
  - Indicate the date of the patient's last menstrual period on the laboratory slip.
- PT** Explain the possibility that serial testing is often required.

### Abnormal findings

#### ▲ Increased levels

Multiple pregnancies  
Placental-site trophoblastic tumor  
Intact molar pregnancy  
Diabetes  
Rh incompatibility

#### ▼ Decreased levels

Placental insufficiency  
Toxemia  
Preeclampsia  
Hydatidiform mole  
Choriocarcinoma

### notes



**human T-cell lymphotropic virus (HTLV) I/II antibody****Type of test** Blood**Normal findings**

Negative

**Test explanation and related physiology**

Several forms of HTLV, a human retrovirus, affect humans. HTLV-I is associated with adult T-cell leukemia or lymphoma. HTLV-II is associated with adult hairy-cell leukemia and neurologic disorders (e.g., tropical spastic paraparesis).

HTLV transmission is similar to HIV transmission (e.g., body fluid contamination, intravenous drug use, sexual contact, breastfeeding).

Blood and organ donors are routinely tested for the presence of anti-HTLV-I/II antibodies. For accurate diagnosis of HTLV-I/II infection, all initially enzyme immunoassay (EIA)-positive results should be verified by a confirmatory test, such as Western blot or line immunoassay. HTLV-I/II can also be directly detected by real-time amplification of the specific HTLV genomic DNA sequences from the blood of infected patients.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**Abnormal findings**

Acute HTLV infection

Adult T-cell leukemia

Hairy-cell leukemia

Tropical spastic paraparesis

**notes**

## 17-hydroxycorticosteroids (17-OCHS)

**Type of test** Urine (24-hour)

### Normal findings

Adult

Male: 3-10 mg/24 hr or 8.3-27.6  $\mu\text{mol/day}$  (SI units)

Female: 2-8 mg/24 hr or 5.2-22.1  $\mu\text{mol/day}$  (SI units)

Elderly: values slightly lower than that of an adult

Children

< 8 years: < 1.5 mg/24 hr

8-12 years: < 4.5 mg/24 hr

### Test explanation and related physiology

Elevated levels of 17-OCHS are seen in patients with hyperfunctioning of the adrenal gland (Cushing syndrome), whether this condition is caused by a pituitary or adrenal tumor, bilateral adrenal hyperplasia, or ectopic tumors producing adrenocorticotropic hormone (ACTH). Low levels of 17-OCHS are seen in patients who have a hypofunctioning adrenal gland (Addison disease) as a result of destruction of the adrenals (by hemorrhage, infarction, metastatic tumor, or autoimmunity), surgical removal of an adrenal gland, congenital enzyme deficiency, hypopituitarism, or adrenal suppression after prolonged exogenous steroid ingestion.

Testing the urine for this hormone metabolite is only an indirect measure of adrenal function. Urine and plasma levels of cortisol (see p. 290) provide a much more accurate measurement of adrenal function. Because the excretion of cortisol metabolites follows a diurnal variation, a 24-hour collection is necessary.

### Interfering factors

- Emotional and physical stress (e.g., infection) and licorice ingestion may cause increased adrenal activity.
- ✔ Drugs that may cause *increased* 17-OCHS levels include acetazolamide, chloral hydrate, chlorpromazine, colchicine, erythromycin, meprobamate, paraldehyde, quinidine, quinine, and spironolactone.
- ✔ Drugs that may cause *decreased* levels include estrogens, oral contraceptives, phenothiazines, and reserpine.

## Procedure and patient care

- See inside front cover for Routine Urine Testing.
- Note that drugs are usually withheld several days before the urine collection. Check with the physician and laboratory for specific guidelines.
- Assess the patient for signs of stress and report these to the physician.

## Abnormal findings

### ▲ Increased levels

Cushing syndrome  
Ectopic ACTH-producing tumors  
Stress  
Adrenal adenoma or carcinoma  
Hyperthyroidism  
Obesity

### ▼ Decreased levels

Adrenal hyperplasia (adrenogenital syndrome)  
Addison disease  
Adrenal suppression from steroid therapy  
Hypopituitarism  
Hypothyroidism

notes

## 5-hydroxyindoleacetic acid (5-HIAA)

**Type of test** Urine (24-hour)

### Normal findings

2-8 mg/24 hr or 10-40  $\mu\text{mol}/\text{day}$  (SI units)

Female levels lower than male levels

### Test explanation and related physiology

Quantitative analysis of urine levels of 5-HIAA is used to detect and follow the clinical course of patients with carcinoid tumors. Carcinoid tumors are serotonin-secreting tumors that may grow in the appendix, intestine, lung, or any tissue derived from the neuroectoderm. These tumors contain *argentaffin-staining (enteroendocrine)* cells, which produce serotonin and other powerful neurohormones that are metabolized by the liver to 5-HIAA and excreted in the urine. These powerful neurohormones are responsible for the clinical presentation of carcinoid syndrome (bronchospasm, flushing, diarrhea).

### Interfering factors

- Foods as listed below in Procedure and patient care
- ✚ Drugs that may cause *increased* 5-HIAA levels include acetaminophen, acetanilid, acetophenetidin, glyceryl guaiacolate, methocarbamol, and reserpine.
- ✚ Drugs that may cause *decreased* levels include aspirin, chlorpromazine, ethyl alcohol, heparin, imipramine, isoniazid, levodopa, monoamine oxidase inhibitors, methenamine, methylidopa, phenothiazines, promethazine, and tricyclic antidepressants.

### Procedure and patient care

- See inside front cover for Routine Urine Testing for 24-hour collection.
- PT Instruct the patient to refrain from eating foods containing serotonin (e.g., plums, pineapples, bananas, eggplant) for several days before and during testing.
- Keep the specimen on ice or in a refrigerator during the 24-hour collection. A preservative is needed.

## Abnormal findings

### ▲ Increased levels

- Carcinoid tumors
- Noncarcinoid illness
- Cystic fibrosis
- Intestinal malabsorption

### ▼ Decreased levels

- Mental depression
- Migraine headaches

notes

## 21-hydroxylase antibodies

**Type of test** Blood

### Normal findings

< 1 U/mL

### Test explanation and related physiology

Chronic primary adrenal insufficiency (Addison disease) is most commonly caused by the insidious autoimmune destruction of the adrenal cortex and is characterized by the presence of adrenal cortex autoantibodies in the serum. It can occur sporadically or in combination with other autoimmune endocrine diseases. This antibody may precipitate this disease. Measurement of this antibody is used in the investigation of causes of adrenal insufficiency.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Autoimmune adrenal insufficiency

Autoimmune polyglandular syndrome

notes

## **hysterosalpingography** (Uterotubography, Uterosalpingography, Hysterogram)

**Type of test** X-ray with contrast dye

### **Normal findings**

Patent fallopian tubes

No defects in uterine cavity

### **Test explanation and related physiology**

In hysterosalpingography, the uterine cavity and fallopian tubes are visualized radiographically after the injection of contrast material through the cervix. Uterine tumors, intrauterine adhesions, and developmental anomalies can be seen. Tubal obstruction caused by internal scarring, tumor, or kinking also can be detected. A possible therapeutic effect of this test is that passage of dye through the tubes may clear mucous plugs, straighten kinked tubes, or break up adhesions. This test also may be used to document adequacy of surgical tubal ligation.

### **Contraindications**

- Patients with infections of the vagina, cervix, or fallopian tubes, because there is risk of extending the infection
- Patients with suspected pregnancy, because contrast material might induce abortion

### **Potential complications**

- Infection of the endometrium (endometritis)
- Infection of the fallopian tubes (salpingitis)
- Uterine perforation
- Allergic reaction to iodinated dye. This rarely occurs because the dye is not administered intravenously

### **Procedure and patient care**

#### **Before**

**PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.

- Determine pregnancy status of the patient.
- Administer sedatives or antispasmodics, if ordered.

**PT** Tell the patient that no food or fluid restrictions are needed.

#### **During**

- Note the following procedural steps:
  1. After voiding, the patient is placed on the fluoroscopy table in the lithotomy position.

## 532 hysterosalpingography

2. With a speculum in the vagina, contrast material is injected through the cervix. The dye fills the entire upper genital tract (uterus and tubes).
  3. Fluoroscopy is performed, and x-ray images are taken.
- Note that this procedure is performed by a physician in approximately 15 to 30 minutes.
- PT** Tell the patient that she may feel occasional transient menstrual-type cramping and that she may have shoulder pain caused by subphrenic irritation from the dye as it leaks into the peritoneal cavity.

### After

- PT** Inform the patient that a vaginal discharge (sometimes bloody) may be present for 1 to 2 days after the test.
- PT** Instruct the patient on the signs and symptoms of infection (e.g., fever, increased pulse rate, pain).

### Abnormal findings

Uterine tumor (e.g., leiomyoma, cancer) or polyps  
Developmental anomaly (e.g., uterus bicornis) of the uterus  
Intrauterine adhesions or polyps  
Uterine fistula  
Obstruction, kinking, or twisting of the fallopian tubes  
Extrauterine pregnancy  
Tumor of the fallopian tubes

### notes



## hysteroscopy


**Type of test** Endoscopy

### Normal findings

Normal structure and function of the uterus

### Test explanation and related physiology

Hysteroscopy is an endoscopic procedure that provides direct visualization of the uterine cavity by inserting a hysteroscope (a thin, telescope-like instrument) through the vagina and cervix and into the uterus (Figure 25). Hysteroscopy can be used to identify the cause of abnormal uterine bleeding, infertility, and

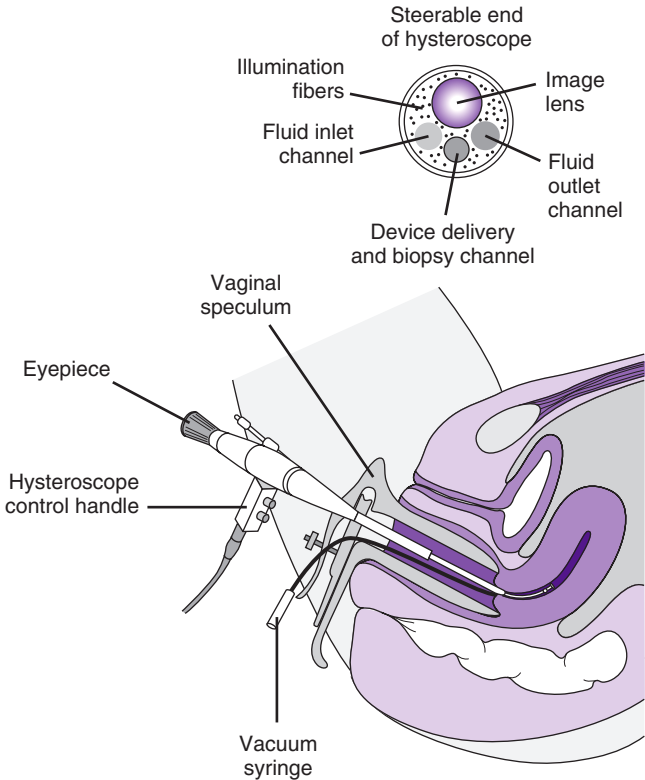

**FIG. 25** Hysteroscopy.

## 534 hysteroscopy

repeated miscarriages. It is also used to identify, evaluate, and perform biopsies of uterine adhesions (Asherman syndrome), polyps, cancer, fibroids, and to detect displaced intrauterine devices (IUDs).

In addition to diagnosing and evaluating uterine problems, hysteroscopy can correct uterine problems. For example, uterine adhesions and small fibroids can be removed through the hysteroscope, thus avoiding open abdominal surgery. Hysteroscopy can also be used to perform endometrial ablation, which destroys the uterine lining to treat some cases of heavy dysfunctional uterine bleeding.

### Contraindications

- Patients with pelvic inflammatory disease
- Patients with vaginal discharge

### Potential complications

- Uterine perforation
- Infection

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Obtain informed consent for this procedure.
- Assess the pregnancy status of the patient.

**PT** Instruct the patient to be NPO (nothing by mouth) for at least 8 hours before the test.

#### During

- Note the following procedural steps:
  1. Hysteroscopy may be performed in the operating room or the doctor's office. Local, regional, general, or no anesthesia may be used. (The type of anesthesia depends on other procedures that may be done at the same time.)
  2. The patient is placed in the lithotomy position. The vaginal area is cleansed with an antiseptic solution.
  3. The cervix may be dilated before this procedure.
  4. The hysteroscope is inserted through the vagina and cervix and into the uterus.
  5. A liquid or gas is released through the hysteroscope to expand the uterus for better visualization.
  6. If minor surgery is to be performed, small instruments will be inserted through the hysteroscope.

7. For more detailed or complicated procedures, a laparoscope may be used to concurrently view the outside of the uterus.
- Note that hysteroscopy is performed by a physician in approximately 30 minutes.

### After

- PT** Tell the patient that it is normal to have slight vaginal bleeding and cramps for a day or two after the procedure.
- PT** Inform the patient that signs of fever, severe abdominal pain, or heavy vaginal discharge or bleeding should be reported to her physician.

### Abnormal findings

Endometrial cancer, polyps, or hyperplasia  
Uterine fibroids  
Asherman syndrome  
Septate uterus  
Displaced IUD

### notes

## immunoglobulin quantification

**Type of test** Blood

### Normal findings

Results vary by age and methods.

#### IgG (mg/dL)

Adults: 565-1765

Children: 250-1600

#### IgA (mg/dL)

Adults: 85-385

Children: 1-350

#### IgM (mg/dL)

Adults: 55-375

Children: 20-200

#### IgD and IgE

Minimal

### Test explanation and related physiology

Antibodies are made up of gamma globulin protein and are called *immunoglobulins*. There are many classes of immunoglobulins. *Immunoglobulin G (IgG)* constitutes approximately 75% of the serum immunoglobulins; therefore it constitutes the majority of circulating blood antibodies. *IgA* constitutes approximately 15% of the immunoglobulins in the body and is present primarily in secretions of the respiratory and gastrointestinal tracts, in saliva, in colostrum, and in tears. IgA is also present to a smaller degree in the blood. *IgM* is primarily responsible for ABO blood grouping and rheumatoid factor; it is also involved in the immunologic reaction to many infections. IgM does not cross the placenta, so an elevation of IgM in a newborn indicates in utero infection (e.g., rubella, cytomegalovirus, or sexually transmitted disease). *IgE* often mediates an allergic response and is measured to detect allergic diseases. *IgD*, which constitutes the smallest part of the immunoglobulins, is rarely evaluated or detected.

Serum protein quantification is used to detect and monitor the course of hypersensitivity diseases, immune deficiencies, autoimmune diseases, chronic infections, and intrauterine fetal infections. Although *electrophoresis* is usually required to interpret an elevated immunoglobulin class as polyclonal versus monoclonal, *immunofixation* is usually required to characterize a monoclonal protein.

Increased serum immunoglobulin concentrations occur because of polyclonal or oligoclonal immunoglobulin proliferation in hepatic disease, connective tissue diseases, and acute and chronic infections. Elevation of immunoglobulins may occur in monoclonal gammopathies (e.g., multiple myeloma, primary systemic amyloidosis, and monoclonal gammopathies of undetermined significance). Decreased immunoglobulin levels are found in patients with acquired or congenital immune deficiencies. It can be used to monitor therapy and recurrence. Testing can determine the type of connective tissue disease, its severity, its clinical course, and its response to therapy.

### Interfering factors

- ✦ Drugs that may cause *increased* immunoglobulin levels include hydralazine, isoniazid, phenytoin, procainamide, tetanus toxoid/antitoxin, and gamma globulin.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Indicate on the laboratory slip if the patient has received any vaccinations or immunizations in the past 6 months.

### Abnormal findings

#### ▲ Increased IgA levels

Chronic liver diseases  
(e.g., primary biliary  
cirrhosis)  
Chronic infections  
Inflammatory bowel disease

#### ▼ Decreased IgA levels

Ataxia or telangiectasia  
Congenital isolated  
deficiency  
Hypoproteinemia (e.g.,  
nephrotic syndrome  
or protein-losing  
enteropathies)  
Immunosuppressive drugs  
(e.g., steroids, dextran)

▲ **Increased IgG levels**

Chronic granulomatous infections (e.g., tuberculosis, Wegener granulomatosis, sarcoidosis)  
Hyperimmunization reactions  
Chronic liver disease  
Multiple myeloma (monoclonal IgG type)  
Autoimmune diseases (e.g., rheumatoid arthritis, Sjögren disease, systemic lupus erythematosus)  
Intrauterine devices

▲ **Increased IgM levels**

Waldenström macroglobulinemia  
Chronic infections (e.g., hepatitis, mononucleosis, sarcoidosis)  
Autoimmune diseases (e.g., systemic lupus erythematosus, rheumatoid arthritis)  
Acute infections  
Chronic liver disorders (e.g., biliary cirrhosis)

▲ **Increased IgE levels**

Allergy reactions (e.g., hay fever, asthma, eczema, anaphylaxis)  
Allergic infections (e.g., aspergillosis or parasites)

▼ **Decreased IgG levels**

Wiskott-Aldrich syndrome  
Agammaglobulinemia  
AIDS  
Hypoproteinemia (e.g., nephrotic syndrome, protein-losing enteropathies)  
Drug immunosuppression (e.g., steroids, dextran)  
Non-IgG multiple myeloma  
Leukemia

▼ **Decreased IgM levels**

Agammaglobulinemia  
AIDS  
Hypoproteinemia (e.g., nephrotic syndrome, protein-losing enteropathies)  
Drug immunosuppression (e.g., steroids, dextran)  
IgG or IgA multiple myeloma  
Leukemia

▼ **Decreased IgE levels**

Agammaglobulinemia

## insulin assay

**Type of test** Blood

### Normal findings

6-26  $\mu\text{U}/\text{mL}$  or 43-186  $\text{pmol}/\text{L}$  (SI units)

Newborn: 3-20  $\mu\text{U}/\text{mL}$

### Possible critical values

> 30  $\mu\text{U}/\text{mL}$

### Test explanation and related physiology

Insulin assay is used to diagnose insulinoma (tumor of the islets of Langerhans). It is also used in the evaluation of patients with fasting hypoglycemia.

Some investigators believe that measuring the ratio of the blood sugar and insulin on the same specimen obtained during the oral glucose tolerance test (GTT; see p. 467) is more reliable than measuring insulin levels alone. Combined with the oral GTT, the insulin assay can show characteristic curves. For example, patients with juvenile diabetes have low fasting insulin levels and display flat GTT insulin curves because of little or no increase in insulin levels. Patients who have mild cases of diabetes have normal fasting insulin levels and display GTT curves with a delayed rise.

Type 2 diabetes (adult onset) is characterized by an excess of insulin production in response to GTT. This hyperresponse of insulin may precede hyperglycemia by many years, allowing the patient time and opportunity to take action to reduce the incidence of outright diabetes through diet management and lifestyle changes.

When combined with fasting blood sugar, the insulin assay is very accurate in detecting insulinoma. After the patient fasts 12 to 14 hours, the insulin/glucose ratio should be less than 0.3. Patients with insulinoma have ratios greater than this. To increase the sensitivity and specificity of these combined tests for insulinoma, Turner and others have proposed the *amended* insulin/glucose ratios using variable mathematical *fudge factors*:

$$\frac{\text{Serum insulin level} \times 100}{\text{Serum glucose} - 30 \text{mg} / 100 \text{mL}}$$

A Turner-amended ratio of more than 50 suggests insulinoma.

### Interfering factors

- Antiinsulin antibodies can interfere with radioimmunoassay.
- Food intake and obesity may cause increased insulin levels.

## 540 insulin assay

- Recent administration of radioisotopes may affect test results.
- ▼ Drugs that may cause *increased* insulin levels include corticosteroids, levodopa, and oral contraceptives.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red
- If the serum insulin level will be measured during the GTT, collect the blood sample before oral ingestion of the glucose load and often at designated intervals after glucose ingestion.

### Abnormal findings

#### ▲ Increased levels

Insulinoma  
Cushing syndrome  
Acromegaly  
Obesity  
Fructose or galactose  
intolerance

#### ▼ Decreased levels

Diabetes mellitus  
Hypopituitarism

### notes



## insulin-like growth factor (IGF-1, Somatomedin C, Insulin-like growth factor binding proteins [IGF BP])

**Type of test** Blood

### Normal findings

Adults: 42-110 ng/mL

Children:

Age (yr)	Girls (ng/mL)	Boys (ng/mL)
0-8	5-128	2-118
9-10	24-158	15-148
11-13	65-226	55-216
14-15	124-242	114-232
16-17	94-231	84-211
18-19	66-186	56-177

### Test explanation and related physiology

Growth hormone (GH) exerts its effects on many tissues through a group of peptides called *somatomedins*. The most commonly tested somatomedins are insulin-like growth factor 1 (IGF-1) and IGF-3. Measurement of free IGF-1 and IGF binding protein (BP) 3 is preferred to GH measurements in cases of short stature in early adolescence. IGF is the test of choice in identifying and monitoring treatment of acromegaly.

Great variation in GH secretion occurs during the day. A random GH assay result may significantly overlap between normal and abnormal values. To diminish the common variations in GH secretion, screening for IGF-1 provides a more accurate reflection of the mean plasma concentration of GH. Somatomedins are not affected (as GH is) by the time of day, food intake, or exercise because they circulate bound to proteins that are durable or long-lasting. Normally there is a large increase during the pubertal growth spurt.

Levels of IGF-1 depend on levels of GH. As a result, IGF-1 levels are low when GH levels are deficient. (See GH [p. 474] for a discussion of causes and diseases associated with GH deficiency.) Nonpituitary causes of reduced IGF-1 levels include malnutrition, severe chronic illnesses, severe liver disease, hypothyroidism, renal failure, inflammatory bowel disease, and Laron dwarfism. Abnormally low test results require an abnormally reduced or absent GH during a GH-stimulation test (see p. 477) to make the diagnosis of GH deficiency.

## 542 insulin-like growth factor

Pediatricians commonly use *IGF BPs* to even further diminish the influence of the variables affecting GH and somatomedin levels. Specifically, IGF BP 2 and IGF BP 3 are the most commonly measured. However, if GH deficiency is strongly suspected yet documentation using GH or somatomedins is questionable, IGF BP determinations are helpful. IGF BP 3 is less age-dependent and is the most accurate (97% sensitivity and specificity). These proteins help to evaluate GH deficiencies and GH-resistant syndromes (e.g., Laron dwarfism). Finally, these binding proteins are very useful in predicting responses to therapeutic exogenous GH administration.

### Interfering factors

- ▼ Drugs that may cause *decreased* levels include high doses of estrogens.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: lavender or red

### Abnormal findings

#### ▲ Increased levels

- Acromegaly
- Gigantism
- Hyperpituitarism
- Obesity
- Pregnancy
- Precocious puberty

#### ▼ Decreased levels

- GH deficiency/resistance
- Laron dwarfism
- Inactive GH
- Resistance to somatomedins
- Nutritional deficiency
- Delayed puberty
- Pituitary tumor
- Hypopituitarism
- Cirrhosis of the liver

notes

## intravascular ultrasound (IVUS)

**Type of test** Ultrasound

### Normal findings

Normal coronary arteries

### Test explanation and related physiology

Percutaneous IVUS imaging requires very small, specially made transducers that are mounted on the tip of an intravascular catheter. Unlike arteriography, which shows a shadow of the arterial lumen, IVUS shows a tomographic, cross-sectional view of the vessel. This orientation enables direct measurements of lumen dimensions, which are considered to be more accurate than angiographic dimensions. The data obtained can be restructured into a longitudinal image by the ultrasound machine software to create a three-dimensional image of the particular segment of artery.

IVUS is an important technology for studying the progression, stabilization, and potential regression of coronary atherosclerosis. IVUS permits imaging of the lumen size, vessel wall structure, and any atheroma that may be present. It allows characterization of atheroma size, plaque distribution, and lesion composition and enables accurate visualization of not only the lumen of the coronary arteries, but also the atheroma that may be hidden within the vessel wall. This procedure is predominantly used in the coronary arteries.

IVUS is used in the following clinical situations:

- Assessment of coronary stent placement and determination of minimum luminal diameter within the stent
- Preinterventional assessment of lesion characteristics and vessel dimensions
- Determination of the mechanism of stent restenosis and selection of appropriate therapy (plaque ablation versus repeat balloon expansion)
- Evaluation of coronary obstruction at a location difficult to image by angiography
- Assessment of a suboptimal angiographic result after stent placement in cases in which the degree of stenosis of a coronary artery is unclear
- Guidance and assessment for vascular atherectomy
- Determination of plaque location and circumference for coronary atherectomy

## 544 intravascular ultrasound

- Determination of the extent of atherosclerosis in patients with characteristic anginal symptoms and a positive functional study with no focal stenoses or mild coronary artery disease on angiography
- Assessment of the changes in plaque volume after lipid-lowering therapy

### Interfering factors

- The accuracy of ultrasonography is dependent on the skills of the sonographer.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that fasting is required.

#### During

- The IVUS probe is placed by coronary angiographic procedures (see p. 201).
- The test is completed in approximately 1 hour, usually by a cardiologist.

#### After

- See cardiac catheterization (see p. 201) for postprocedure care.

### Abnormal findings

Coronary occlusive disease

notes

## intrinsic factor antibody (IF ab)

**Type of test** Blood

### Normal findings

Negative

### Test explanation and related physiology

IF ab is used to diagnose pernicious anemia (PA). PA is one of the major causes of vitamin B<sub>12</sub> deficiency and megaloblastic anemia. It is a disease of the stomach in which secretion of IF is severely reduced or absent, resulting in malabsorption of vitamin B<sub>12</sub>. PA appears to be an autoimmune process.

Approximately 50% to 75% of adult patients have IF abs. There are two types of this antibody. Type I, blocking antibody, is more common and prevents the binding of vitamin B<sub>12</sub> and IF. Type II binding antibody is less specific for PA and affects the binding of IF in the ileum. The blocking antibody is extremely specific for PA and is more sensitive than the binding antibody. In the context of a low or borderline B<sub>12</sub> result, in which other clinical and hematologic findings are compatible with a diagnosis of B<sub>12</sub> deficiency, the presence of IF blocking antibody can be taken as confirmation of this diagnosis and, at the same time, as an indication of its cause. A negative result, on the other hand, cannot rule out the possibility of PA because blocking antibody is not demonstrable in nearly 50% of all patients with this disorder.

### Interfering factors

- IF ab levels are *decreased* if an injection of vitamin B<sub>12</sub> is administered within 48 hours of testing.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Pernicious anemia

notes

**iron level and total iron-binding capacity** (Fe and TIBC, Transferrin saturation, Transferrin)**Type of test** Blood**Normal findings****Iron**Male: 80-180 mcg/dL or 14-32  $\mu\text{mol/L}$  (SI units)Female: 60-160 mcg/dL or 11-29  $\mu\text{mol/L}$  (SI units)

Newborn: 100-250 mcg/dL

Child: 50-120 mcg/dL

**TIBC**250-460 mcg/dL or 45-82  $\mu\text{mol/L}$  (SI units)**Transferrin**

Adult male: 215-365 mg/dL or 2.15-3.65 g/L (SI units)

Adult female: 250-380 mg/dL or 2.50-3.80 g/L (SI units)

Newborn: 130-275 mg/dL

Child: 203-360 mg/dL

**Transferrin saturation**

Male: 20%-50%

Female: 15%-50%

**Test explanation and related physiology****Serum iron**

Abnormal levels of iron are characteristic of many diseases, including iron deficiency anemia and hemochromatosis. Seventy percent of iron in the body is found in the hemoglobin of red blood cells (RBCs). The other 30% is stored iron in the form of *ferritin* (see p. 404) and hemosiderin. Iron is supplied by the diet. Iron is bound to a globulin protein called *transferrin*. When iron stores are low, transferrin levels increase. Transferrin is low when there is too much iron. Usually about one-third of the transferrin is used to transport iron. Because of this, the blood serum has considerable extra iron-binding capacity, which is the *unsaturated iron-binding capacity (UIBC)*. The TIBC equals UIBC plus the serum iron measurement. Some laboratories measure UIBC, some measure TIBC, and some measure transferrin. The serum iron determination is a measurement of the quantity of iron bound to transferrin.

Iron-deficiency anemia is a result of reduced serum iron. Iron-deficiency anemia has many causes, including the following:

- Insufficient iron intake
- Inadequate gut absorption
- Increased requirements (e.g., in growing children)
- Loss of blood (e.g., menstruation, bleeding peptic ulcer)

Iron deficiency results in decreased production of hemoglobin, which in turn results in small, pale (microcytic, hypochromic) RBCs. A decreased serum iron level, elevated TIBC, and low transferrin saturation (TS) value are characteristic of iron deficiency anemia.

Acute iron poisoning due to accidental or intentional overdose is characterized by a serum iron level that exceeds the TIBC. Chronic iron overload or poisoning is called *hemochromatosis* or *hemosiderosis*. Excess iron is usually deposited in the brain, liver, and heart and causes severe dysfunction of these organs. Massive blood transfusions also may cause elevated serum iron levels, although only transiently.

### Total iron-binding capacity and transferrin

TIBC is a measurement of all proteins available for binding mobile iron. Transferrin represents the largest quantity of iron-binding proteins. Therefore TIBC is an indirect yet accurate measurement of transferrin. Ferritin is not included in TIBC because it binds only stored iron. TIBC is increased in 70% of patients with iron deficiency.

Transferrin is a *negative* acute-phase reactant protein. That is, in various acute inflammatory reactions, transferrin levels diminish. Transferrin also is diminished in the face of chronic illnesses, such as malignancy, collagen vascular diseases, or liver diseases. Hypoproteinemia is also associated with reduced transferrin levels.

TIBC varies minimally according to iron intake and is more of a reflection of liver function (transferrin is produced by the liver) and nutrition than of iron metabolism.

### Total iron-binding capacity and transferrin saturation

The percentage of transferrin and other mobile iron-binding proteins saturated with iron is calculated by dividing the serum iron level by the TIBC:

$$TS(\%) = \frac{\text{Serum iron level}}{\text{TIBC}} \times 100\%$$

The normal value for TS is 20% to 50%. TS is decreased to less than 15% in patients with iron-deficiency anemia. It is increased in

## 548 iron level and total iron-binding capacity

patients with hemolytic, sideroblastic, or megaloblastic anemias. TS is also increased in patients with iron overload or poisoning. Increased intake or absorption of iron (as in hemochromatosis) leads to elevated iron levels. In such cases, TIBC is unchanged; as a result, the percentage of TS is very high. UIBC has been proposed as an inexpensive alternative to transferrin saturation.

Chronic illness is characterized by a low serum iron level, decreased TIBC, and normal TS. Pregnancy is marked by high levels of protein, including transferrin. Because iron requirements are high, it is not unusual to find low serum iron levels, high TIBC, and a low percentage of TS in late pregnancy.

### Contraindications

- Patients with hemolytic diseases, because they may have an artificially high iron content

### Interfering factors

- Recent blood transfusions may affect test results.
- Recent ingestion of a meal containing high iron content may affect test results.
- Hemolytic diseases may be associated with an artificially high iron content.
- ✚ Drugs that may cause *increased* iron levels include chloramphenicol, dextran, estrogens, ethanol, iron preparations, methyldopa, and oral contraceptives.
- ✚ Drugs that may cause *decreased* iron levels include adrenocorticotrophic hormone (ACTH), chloramphenicol, cholestyramine, colchicine, deferoxamine, methicillin, and testosterone.
- ✚ Drugs that may cause *increased* TIBC levels include fluorides and oral contraceptives.
- ✚ Drugs that may cause *decreased* TIBC levels include ACTH and chloramphenicol.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes (12 hours)
- Blood tube commonly used: red
- Avoid hemolysis, because the iron contained in the RBCs will pour out into the serum and cause artificially high iron levels.



**Abnormal findings****▲ Increased serum iron levels**

Hemosiderosis  
 Hemochromatosis  
 Hemolytic anemia  
 Hepatitis  
 Hepatic necrosis  
 Lead toxicity  
 Iron poisoning  
 Massive blood transfusion

**▲ Increased TIBC or transferrin levels**

Oral contraceptives  
 Pregnancy (late)  
 Polycythemia vera  
 Iron deficiency anemia

**▲ Increased TS or TIBC saturation**

Hemochromatosis  
 Hemosiderosis  
 Acute iron overdose  
 Hemolytic anemia

**▼ Decreased serum iron levels**

Insufficient dietary iron  
 Chronic blood loss  
 Inadequate absorption of iron  
 Pregnancy (late)  
 Iron deficiency anemia  
 Neoplasia  
 Chronic gastrointestinal blood loss  
 Chronic hematuria  
 Chronic heavy physiologic or pathologic menstruation

**▼ Decreased TIBC or transferrin levels**

Hypoproteinemia  
 Inflammatory diseases  
 Cirrhosis  
 Hemolytic anemia  
 Pernicious anemia  
 Sickle cell anemia

**▼ Decreased TS or TIBC saturation**

Iron deficiency anemia  
 Chronic illnesses (e.g., malignancy)

notes

**ischemia-modified albumin (IMA)****Type of test** Blood**Normal findings**

&lt; 85 IU/mL

**Test explanation and related physiology**

When albumin is exposed to an ischemic environment, its N terminal is altered; this causes an alteration of the albumin called *ischemia-modified albumin (IMA)*. This has become particularly helpful in identifying cardiac ischemia in patients with chest pain. When combined with troponins (see p. 912), myoglobin (see p. 635), and electrocardiography, the diagnosis of an ischemic cardiac event can be corroborated or ruled out. IMA is produced continually during the period of ischemia. Blood levels rise within 10 minutes of the initiation of the ischemic event and stay elevated for 6 hours after ischemia has resolved.

IMA may also be elevated in patients with pulmonary embolus or acute stroke. False positives can occur in other clinical circumstances, such as advanced cancers, acute infections, and end-stage renal or liver disease.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: yellow
- This test is usually done after the initial onset of chest pain, then 12 hours later, and then daily testing for 3 to 5 days.
- Record the exact time and date of venipuncture on each laboratory slip. This aids in the interpretation of the temporal pattern of blood level elevations.

**Abnormal findings****▲ Increased levels**

Myocardial ischemia  
Brain ischemia  
Pulmonary ischemia

notes

## lactic acid (Lactate)

**Type of test** Blood

### Normal findings

Venous blood: 5-20 mg/dL or 0.6-2.2 mmol/L (SI units)

Arterial blood: 3-7 mg/dL or 0.3-0.8 mmol/L (SI units)

### Test explanation and related physiology

Under conditions of normal oxygen availability to tissues, glucose is metabolized to  $\text{CO}_2$  and  $\text{H}_2\text{O}$  for energy. When oxygen to the tissues is diminished, anaerobic metabolism of glucose occurs, and lactate (lactic acid) is formed instead of  $\text{CO}_2$  and  $\text{H}_2\text{O}$ . To compound the problem of lactic acid buildup, when the liver is hypoxic, it fails to clear the lactic acid. Lactic acid levels accumulate, causing lactic acidosis (LA). Lactic acid blood levels are used to document the presence of tissue hypoxia, determine the degree of hypoxia, and monitor the effect of therapy.

Type I LA is caused by diseases or factors that increase lactate but are not hypoxic related (e.g., liver diseases or drugs). LA caused by hypoxia is classified as type II. Shock, convulsions, and extremity ischemia are the most common causes of type II LA. Type III LA is idiopathic and is most commonly seen in nonketotic patients with diabetes.

### Interfering factors

- The prolonged use of a tourniquet or clenching of hands increases lactate levels.
- ✚ Drugs that *increase* levels include aspirin, cyanide, ethanol (chronic use), nalidixic acid, and phenformin.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: no
  - Blood tube commonly used: red
- PT** Instruct the patient to avoid making a fist before and while blood is being withdrawn.
- Avoid the use of a tourniquet if possible.

## Abnormal findings

### ▲ Increased levels

Shock

Tissue ischemia

Carbon monoxide poisoning

Severe liver disease

Genetic errors of metabolism

Diabetes mellitus (nonketotic)

notes

## lactic dehydrogenase (LDH, Lactate dehydrogenase)

**Type of test** Blood

### Normal findings

#### Total

Adult/elderly: 100-190 units/L at 37° C (lactate → pyruvate) or  
100-190 units/L (SI units)

Child: 60-170 units/L (30° C)

Infant: 100-250 units/L

Newborn: 160-450 units/L

#### Isoenzymes

Adult/elderly:

LDH-1: 17%-27%

LDH-2: 27%-37%

LDH-3: 18%-25%

LDH-4: 3%-8%

LDH-5: 0%-5%

### Test explanation and related physiology

The LDH test is a measure of total LDH. There are actually five separate fractions (isoenzymes) that make up the total LDH. In general, isoenzyme LDH-1 comes mainly from the heart; LDH-2 comes primarily from the reticuloendothelial system; LDH-3 comes from the lungs and other tissues; LDH-4 comes from the kidney, placenta, and pancreas; and LDH-5 comes mainly from the liver and striated muscle. In normal persons, LDH-2 makes up the greatest percentage of total LDH.

With myocardial injury, the serum LDH level rises within 24 to 48 hours after a myocardial infarction (MI), peaks in 2 to 3 days, and returns to normal in approximately 5 to 10 days. Other cardiac markers (e.g., CK-MB, p. 297; and troponin, p. 912) have replaced the indications for LDH in patients with MI. LDH-1 is generally not as useful as troponin or creatine kinase-MB for detection of MI unless the MI occurred 24 hours or more before the assay. LDH-5 is usually not as reliable as the transaminases as a liver function test.

LDH is also measured in other body fluids. Elevated urine levels of total LDH indicate neoplasm or injury to the urologic system. When the LDH in an effusion (pleural, cardiac, or peritoneal) is more than 60% of the serum total LDH (i.e., effusion LDH/serum LDH ratio is > 0.6), the effusion is said to be an *exudate* and not a *transudate*.

## 554 lactic dehydrogenase

### Interfering factors

- Strenuous exercise may cause an elevation of total LDH, specifically LDH-1, -2, and -5.
- Hemolysis of blood will cause false-positive LDH levels.
- ✚ Drugs that may cause *increased* LDH levels include alcohol, anesthetics, aspirin, clofibrate, fluorides, mithramycin, narcotics, and procainamide.
- ✚ Drugs that may cause *decreased* levels include ascorbic acid.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased values

Myocardial infarction

Pulmonary disease (e.g., embolism, infarction, pneumonia)

Hepatic disease (e.g., hepatitis, active cirrhosis, neoplasm)

Red blood cell disease (e.g., hemolytic or megaloblastic anemia or red blood cell destruction from prosthetic heart valves)

Skeletal muscle disease and injury (e.g., muscular trauma)

Renal parenchymal disease (e.g., infarction, glomerulonephritis, acute tubular necrosis)

Intestinal ischemia and infarction

Testicular tumors (seminoma or dysgerminomas)

Lymphoma and other reticuloendothelial system tumors

Advanced solid tumor malignancies

Pancreatitis

Diffuse disease or injury (e.g., heatstroke)

### notes

## lactoferrin

**Type of test** Stool

### Normal findings

None detected

### Test explanation and related physiology

Lactoferrin is a glycoprotein expressed by activated neutrophils. The detection of lactoferrin in a fecal sample therefore serves as a surrogate marker for inflammatory white blood cells (WBCs) in the intestinal tract. WBCs in the stool are not stable and may be easily destroyed by temperature changes, delays in testing, and toxins within the stool. As a result, WBCs may not be detected by common microscopic methods. Lactoferrin assay has allowed the identification of inflammatory cells in the stool without the use of microscopy.

Detection of fecal lactoferrin allows for the differentiation of inflammatory and noninflammatory intestinal disorders in patients with diarrhea. Usually the test is used as a diagnostic aid to help identify patients with active inflammatory bowel disease (e.g., Crohn disease or ulcerative colitis) and to rule out those with active irritable bowel syndrome, which is non-inflammatory. Lactoferrin is also present in patients with bacterial enteritis (e.g., *Shigella* spp., *Salmonella* spp., *Campylobacter jejuni*, and *Clostridium difficile*). Diarrhea caused by viruses and most parasites is not associated with elevated lactoferrin levels.

### Interfering factors

- The stool specimen should be examined immediately. In some instances, a specific stool preservative–enteric transport media can be used.
- Breastfeeding can affect test results in breastfed infants.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

**PT** Instruct the patient not to mix urine or toilet paper with the specimen.

#### During

- Stool is collected in a clean bedpan.
- Place at least 5 g of stool in a clean specimen container.

## 556 lactoferrin

### After

- Observe appropriate contamination precautions.
  - Transfer the specimen to the laboratory immediately.
- PT** Inform the patient that results are usually available in less than a half hour.

### Abnormal findings

Bacterial enteritis

Crohn disease

Ulcerative colitis

### notes



## lactose tolerance test (Hydrogen breath test)

**Type of test** Blood

### Normal findings

**Blood:** Adult/elderly: rise in plasma glucose levels  $> 20$  mg/dL  
No abdominal cramps or diarrhea

**Breath:**  $< 50$  ppm hydrogen increase over baseline

### Test explanation and related physiology

This test is performed to detect lactose intolerance, intestinal malabsorption, maldigestion, or bacterial overgrowth in the small intestine. Because lactose-intolerant patients have an absence of lactase, any lactose (the common sugar in dairy products) ingested will not be digested in the small bowel. Thus the colon is flooded with a high lactose load. Although all adults have some degree of lactase reduction, severe lactose intolerance can occur in patients with inflammatory bowel disease, short-gut syndrome, and other malabsorption syndromes. Lactase deficiency can be congenital and become apparent in newborns. These infants present with vomiting, diarrhea, malabsorption, and failure to thrive.

In this test, the patient is given an oral lactose load. If lactase is not present in sufficient quantities, lactose is not metabolized to glucose and galactose. Plasma levels of glucose do not rise as expected. Therefore lower than expected serum glucose levels suggest intestinal lactase deficiency. Patients who have malabsorption without lactase deficiency also fail to elevate the blood glucose levels.

There is a breath-test portion of this diagnostic test in which exhaled air is analyzed for hydrogen content. This is called the *lactose breath test* (or *hydrogen breath test*). The bacteria in the colon produce hydrogen when exposed to unabsorbed food, particularly the lactose load that was not absorbed in the small intestine. Large amounts of hydrogen may also be produced when the colonic bacteria move back into the small intestine, a condition called *bacterial overgrowth of the small bowel*. In this instance the overgrowth bacteria are exposed to the lactose load that has not had a chance to completely traverse the small intestine to be fully digested and absorbed. Large amounts of the hydrogen produced by the bacteria are absorbed into the blood flowing through the wall of the small intestine and colon. This hydrogen-containing blood travels to the lungs, where the hydrogen is released and exhaled in the breath in measurable quantities.

## 558 lactose tolerance test

Before lactose hydrogen breath testing, individuals must fast for at least 12 hours. At the start of the test, the individual blows into a hydrogen analyzer. The individual then ingests a small amount of the test sugar (e.g., lactose, sucrose, sorbitol, fructose, or lactulose, depending on the purpose of the test). Additional samples of breath are collected and analyzed for hydrogen every 15 minutes for 1 to 5 hours. When rapid intestinal transit is present, the test dose of nondigestible lactulose reaches the colon more quickly than normal; therefore hydrogen is produced by the colonic bacteria soon after the sugar is ingested. When bacterial overgrowth of the small bowel is present, ingestion of lactulose results in two separate periods during the test in which hydrogen is produced: an earlier period caused by the bacteria in the small intestine and a later one caused by the bacteria in the colon.

### Interfering factors

- Enterogenous steatorrhea
- Strenuous exercise
- Smoking may increase blood glucose levels.
- Ethnicity has a major effect on primary lactose deficiency.
- Patients with diabetes may have glucose levels that exceed 20 mg/dL despite lactase insufficiency.
- ✚ Antibiotics can *decrease* the bacteria in the intestine and may cause false-negative breath tests and thus should not be taken for 1 month before testing.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. Inform the patient that four blood samples will be needed.
- PT** Instruct the patient to fast for 12 hours before testing.
- PT** Instruct the patient to avoid strenuous exercise for 8 hours before testing because exercise may factitiously affect the blood glucose level.
- PT** Inform the patient that smoking is prohibited for approximately 8 hours before testing because smoking can increase the blood glucose level.

#### During

- Collect a venous blood sample in a gray-top tube from the fasting patient.
- Provide a specified dose of lactose for the patient. Usually dilute 50 to 100 g of lactose with 200 mL of water for ingestion in adults.
- Note that pediatric doses of lactose are based on weight.

- Collect three more blood samples at 30, 60, and 120 minutes after the ingestion of lactose.
- PT Tell the patient that the only discomfort is the venipuncture; however, patients with lactase deficiency may have symptoms of lactose intolerance (e.g., cramps and diarrhea).
- If the breath test is being done, the exhaled air is evaluated for hydrogen content before ingestion of lactose and every 15 minutes thereafter. Hydrogen levels are recorded.

### After

- Apply pressure to the venipuncture site.
- Note that patients with abnormal test results may receive a monosaccharide tolerance test (e.g., glucose or galactose tolerance test).

### Abnormal findings

#### ▼ Decreased levels

Lactase insufficiency

Intestinal malabsorption or maldigestion

Small bowel overgrowth of bacteria

notes

## Legionnaires disease antibody test

**Type of test** Blood

### Normal findings

No *Legionella* antibody titer

### Test explanation and related physiology

Legionnaires disease was originally described as a fulminating pneumonia caused by *Legionella pneumophila*, a tiny gram-negative, rod-shaped bacterium. This organism can also cause an influenza-type illness called *Pontiac fever*.

The diagnosis of Legionnaires disease can be made by culturing this organism from suspected infected fluids (e.g., blood, sputum, lung tissue, and pleural fluid). Sputum for this test is best obtained by transtracheal aspiration or from bronchial washings. Another method of diagnosis is by directly identifying the organism in a microscopic smear of infected fluid with the use of direct fluorescent antibody methods. If positive, this allows for a rapid diagnosis of *Legionella*.

The most common and easiest method for diagnosis is detection of the antibody directed against the *Legionella* bacterium in the patient's blood. A presumptive diagnosis of Legionnaires disease can be made in a symptomatic person when a single antibody titer is 1:256 or greater. A fourfold rise in titer to at least 1:128 between the acute-phase (1-week) and the convalescent-phase (3-week) titer is diagnostic.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- For culture: Obtain sputum as indicated in sputum culture (see p. 841).

### Abnormal findings

#### ▲ Increased levels

Legionnaires disease

## leucine aminopeptidase (LAP)

**Type of test** Blood; urine (24-hour)

### Normal findings

#### Blood

Male: 80-200 units/mL or 19.2-48 units/L (SI units)

Female: 75-185 units/mL or 18.0-44.4 units/L (SI units)

#### Urine

2-18 units/24 hr

### Test explanation and related physiology

LAP is an intracellular enzyme that exists in the hepatobiliary system and, to a much smaller degree, in the pancreas and small intestine. LAP is mainly used in diagnosing liver disorders and in the differential diagnosis of increased levels of alkaline phosphatase (ALP; see p. 29). LAP levels tend to parallel ALP levels in hepatic disease. LAP is a sensitive indicator of cholestasis; however, unlike ALP, LAP remains normal in bone disease.

### Interfering factors

- Pregnancy may cause increased values.
- ✦ Estrogens and progesterones may cause *increased* LAP levels.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- If a urine sample is needed, see inside front cover for Routine Urine Testing.

### Abnormal findings

#### ▲ Increased levels

Hepatitis

Cirrhosis

Hepatic necrosis, ischemia, tumor

Hepatotoxic drugs

Cholestasis

Gallstones

### notes

## lipase

**Type of test** Blood

### Normal findings

0-160 units/L or 0-160 units/L (SI units) (Values are method dependent.)

### Test explanation and related physiology

The most common cause of an elevated serum lipase level is acute pancreatitis. Lipase is an enzyme secreted by the pancreas into the duodenum to break down triglycerides into fatty acids. As with amylase (see p. 56), lipase appears in the bloodstream after damage to or disease affecting the pancreatic acinar cells.

Because lipase was thought to be produced only in the pancreas, elevated serum levels were considered to be specific to pathologic pancreatic conditions. It is now apparent that other conditions can be associated with elevated lipase levels. Lipase is excreted through the kidneys. Therefore elevated lipase levels are often found in patients with renal failure. Intestinal infarction or obstruction also can be associated with lipase elevation. However, the lipase elevations in nonpancreatic diseases are less than three times the upper limit of normal compared with those in pancreatitis, in which they are often 5 to 10 times normal values.

In acute pancreatitis, elevated lipase levels usually parallel serum amylase levels. The lipase levels usually rise a little later than amylase (24-48 hours after the onset of pancreatitis) and remain elevated for 5 to 7 days. Because lipase peaks later and remains elevated longer than serum amylase, it is more useful in the diagnosis of acute pancreatitis later in the course of the disease. Lipase levels are less useful in more chronic pancreatic diseases (e.g., chronic pancreatitis, pancreatic carcinoma).

### Interfering factors

- ❖ Drugs that may cause *increased* lipase levels include bethanechol, cholinergics, codeine, indomethacin, meperidine, methacholine, and morphine.
- ❖ Drugs that may cause *decreased* levels include calcium ions.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red

## Abnormal findings

Acute pancreatitis  
Chronic relapsing pancreatitis  
Pancreatic cancer  
Pancreatic pseudocyst  
Acute cholecystitis  
Cholangitis  
Extrahepatic duct obstruction  
Renal failure  
Bowel obstruction or infarction  
Salivary gland inflammation or tumor  
Peptic ulcer disease

## notes

**lipoprotein-associated phospholipase A<sub>2</sub>**  
(Lp-PLA<sub>2</sub> PLAC test)**Type of test** Blood**Normal findings**

Average value for females: 174 ng/mL (range: 120-342)

Average value for males: 251 ng/mL (range: 131-376)

**Test explanation and related physiology**

Lipoprotein-associated phospholipase A<sub>2</sub> (Lp-PLA<sub>2</sub>) promotes vascular inflammation contributing directly to the atherogenic process. Lp-PLA<sub>2</sub> is an independent predictor of cardiovascular disease. When combined with C-reactive protein (CRP) (see p. 295), testing for Lp-PLA<sub>2</sub> markedly increases the predictive value in determining risks for a cardiac event, especially in patients whose Adult Treatment Panel III cardiac risks are moderate. An Lp-PLA<sub>2</sub> level greater than 200 ng/mL would warrant reclassifying the patient to the next highest risk category, which would require more aggressive use of cholesterol-lowering agents. Lp-PLA<sub>2</sub> may play an important role in the progression of atherosclerosis and overall plaque stability.

Lp-PLA<sub>2</sub> is also an accurate aid in assessing the risk for ischemic stroke associated with atherosclerosis at all levels of blood pressure. The PLAC test is an enzyme-linked immunosorbent assay (ELISA) using two highly specific monoclonal antibodies to measure the level of Lp-PLA<sub>2</sub> in the blood.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**Abnormal findings****▲ Increased levels**

Atherosclerosis

**notes**



**lipoproteins** (High-density lipoproteins [HDLs, HDL-C], Low-density lipoproteins [LDLs, LDL-C], Very-low-density lipoproteins [VLDLs], Lipoprotein electrophoresis, Lipoprotein phenotyping, Lipid fractionation, Non-HDL cholesterol, Lipid profile)

**Type of test** Blood

### Normal findings

#### HDL

Male: > 45 mg/dL or > 0.75 mmol/L (SI units)

Female: > 55 mg/dL or > 0.91 mmol/L (SI units)

#### LDL

Adult: < 130 mg/dL

Children: < 110 mg/dL

#### VLDL

7-32 mg/dL

### Test explanation and related physiology

Lipoproteins are considered to be an accurate predictor of heart disease. As part of the lipid profile, these tests are performed to identify persons at risk for developing heart disease and to monitor the response to therapy if abnormalities are found.

Lipoproteins are proteins that are used as markers indicating the levels of lipids within the bloodstream. The lipid profile usually measures total cholesterol, triglycerides, HDL, LDL, and VLDL. HDLs (good cholesterol) are carriers of cholesterol. They are produced in the liver and, to a smaller degree, in the intestines. In general, the purpose of HDLs is believed to be removal of the cholesterol from the peripheral tissues and transportation to the liver for excretion. The function of removing lipids from the endothelium (reverse cholesterol transport) provides a protective effect against heart disease.

Total HDL cholesterol is an independent, inverse risk factor for coronary artery disease (CAD). Whereas low levels (< 35 mg/dL) are believed to increase a person's risk for CAD, high levels (> 60 mg/dL) are considered protective. When HDL and total cholesterol measurements are combined in a ratio fashion, the accuracy of predicting CAD is increased. The total cholesterol/HDL ratio should be at least 5:1, with 3:1 being ideal.

Five subclasses of HDL (2a, 2b, 3a, 3b, and 3c) have been identified, but only 2b is cardioprotective. HDL 2b is the most efficient form of HDL in reverse cholesterol transport. Patients with low total HDL levels often have low levels of HDL 2b. When levels of total HDL are between 40 and 60 mg/dL,

cardioprotective levels of HDL 2b are minimal. However, when levels of total HDL are greater than 60 mg/dL, levels of HDL 2b predominate, and efficient reverse cholesterol transport takes place. This protects the coronary arteries from disease. The other subclasses of HDL are not capable of reverse cholesterol transport and therefore are not cardioprotective.

LDLs (bad cholesterol) are also cholesterol rich. However, most cholesterol carried by LDLs can be deposited into the lining of the blood vessels and is associated with an increased risk of arteriosclerotic heart and peripheral vascular disease. Therefore high levels of LDLs are atherogenic. The LDL level should be less than 70 mg/dL in patients at high risk for heart disease. For patients at moderately high risk, the LDL should be less than or equal to 100 mg/dL (depending on other cardiac risk factors). LDL can be difficult to measure. Therefore the LDL level can be calculated using the Friedewald formula. In this formula, subtracting the HDL plus one-fifth of the triglycerides from the total cholesterol derives LDL level:

$$\text{LDL} = \text{total cholesterol} - (\text{HDL} + [\text{triglycerides} \div 5])$$

There are other formulas for deriving LDL, which may account for different sets of normal values. Furthermore, the formula is inaccurate if the triglycerides exceed 400 mg/dL.

LDL has been divided into seven classes based on particle size. These subclasses are (from largest to smallest) LDL I, LDL IIa, LDL IIb, LDL IIIa, LDL IIIb, LDL IVa, and LDL IVb. The most commonly elevated forms of LDL (IIIa and IIIb) are small enough to get between the endothelial cells and cause atherosclerotic disease. The larger LDL particles (LDL I, LDL IIa, and LDL IIb) cannot get into the endothelial layer and therefore are not associated with increased risk of disease. LDL IVa and IVb, however, are very small and are associated with aggressive arterial plaques that are particularly vulnerable to ulceration and vascular occlusion. Nearly all patients with levels of LDL IVa and IVb greater than 10% of total LDL have vascular events within months.

LDL patterns can be identified, and they are associated with variable risks of CAD. LDL pattern A is seen in patients with mostly large LDL particles and does not carry increased risks for CAD. LDL pattern B is seen in patients with mostly small LDL particles and is associated with an increased risk of CAD. An intermediate pattern is noted in a large number of patients; they have small and large LDL particles and experience an intermediate risk of CAD. LDL levels can be lowered with diet, exercise, and statins.

VLDLs, although carrying a small amount of cholesterol, are the predominant carriers of blood triglycerides. To a lesser degree, VLDLs are also associated with an increased risk of CAD by virtue of their capability to be converted to LDL by lipoprotein lipase in skeletal muscle. The VLDL value is sometimes expressed as a percentage of total cholesterol. Levels in excess of 25% to 50% are associated with increased risk of coronary disease.

### Interfering factors

- Smoking and alcohol ingestion decrease HDL levels.
- Binge eating can alter lipoprotein values.
- HDL values are age and sex dependent.
- HDL values, like cholesterol, tend to decrease significantly for as long as 3 months after MI.
- HDL is elevated in hypothyroid patients and diminished in hyperthyroid patients.
- High triglyceride levels can make LDL calculations inaccurate.
- ▼ Drugs that may cause altered lipoprotein levels include alpha-blockers, aspirin, beta-blockers, estrogens, phenothiazines, phenytoin, steroids, and sulfonamides.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes (12-14 hours)
- Blood tube commonly used: red
- PT Inform the patient that dietary indiscretion within the previous few weeks may influence lipoprotein levels.
- PT Instruct patients with high lipoprotein levels regarding diet, exercise, and appropriate body weight.

### Abnormal findings

#### ▲ Increased HDL levels

Familial HDL  
lipoproteinemia  
Excessive exercise

#### ▼ Decreased HDL levels

Metabolic syndrome  
Familial low HDL  
Hepatocellular disease (e.g.,  
hepatitis or cirrhosis)  
Hypoproteinemia (e.g.,  
nephrotic syndrome or  
malnutrition)

## 568 lipoproteins

### ▲ Increased LDL and VLDL levels

Familial LDL lipoproteinemia  
Nephrotic syndrome  
Glycogen storage diseases (e.g., von Gierke disease)  
Hypothyroidism  
Alcohol consumption  
Chronic liver disease (e.g., hepatitis or cirrhosis)  
Hepatoma  
Gammopathies (e.g., multiple myeloma)  
Familial hypercholesterolemia type IIa  
Cushing syndrome  
Apoprotein CII deficiency

### ▼ Decreased LDL and VLDL levels

Familial hypolipoproteinemia  
Hypoproteinemia (e.g., malabsorption, severe burns, or malnutrition)  
Hyperthyroidism

## notes

**liquid biopsy** (Fluid biopsy, Fluid phase biopsy)

**Type of test** Blood, urine, cerebrospinal fluid (CSF)

**Normal findings**

No abnormal cells or cell nucleic acids

**Test explanation and related physiology**

The principle supporting liquid biopsy is based on repeated findings that tumors shed molecules and cells into bodily fluids (particularly blood). Since the development of molecular testing laboratory methods, researchers have shown that analyzing these molecules and cells can reveal some of the same and even additional information that tissue biopsies provide. With this additional information, personalized treatment strategies can be developed.

Tissue biopsies usually involve an invasive procedure that can be risky, costly, painful, and not easily repeated during the course of the disease therapy. Furthermore, some tumors may be inaccessible. Liquid biopsies are easily obtained through venipuncture or by obtaining a urine specimen. They can track tumors as they develop and change over time.

A liquid biopsy may be used to help find cancer at an early stage. Being able to take multiple samples of blood over time may also help doctors understand what kind of molecular changes are taking place in a tumor. Repeated liquid biopsies can be used in surveillance of treated cancer patients to look for recurrence. In a similar manner, repeated liquid biopsies during treatment can more accurately assess and monitor treatment response.

It is important, at this time, however, to recognize that liquid biopsy is far from being “standard of care.” There are many questions that must be more fully answered before this form of testing becomes clinically commonplace.

Different liquid biopsy tests analyze different kinds of tumor material, such as entire circulating free tumor cells (CTCs), circulating tumor DNA (ctDNA), RNA, proteins, and tiny vesicles of nucleic acid molecules called exomes. The tests detect these molecules or cells in various bodily fluids, including blood, urine, CSF, and saliva.

Liquid biopsies are now being used in lung and breast cancer patients in order to help direct more specific tumor-related therapies. In cardiology, endothelial cells identified in the blood,

## 570 liquid biopsy

may indicate myocardial infarction (MI). Cell free fetal DNA (see cell free maternal testing, p. 217) is used to identify serious fetal genetic and chromosomal abnormalities.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: see laboratory for directions

### Abnormal findings

Presence of circulating tumor cells

Presence of DNA or RNA fragments of known tumor cells

notes

## liver biopsy

**Type of test** Microscopic examination of tissue


### Normal findings

Normal liver histology

### Test explanation and related physiology

Liver biopsy is a safe, simple, and valuable method of diagnosing pathologic liver conditions. For this study, a specially designed needle is inserted through the abdominal wall and into the liver (Figure 26). A piece of liver tissue is removed for microscopic examination. Percutaneous liver biopsy is used in the diagnosis of various liver disorders (e.g., cirrhosis, hepatitis, drug reaction, granuloma, and tumor). Biopsy is indicated for the following:

- Patients with unexplained hepatomegaly
- Patients with persistently elevated liver enzymes

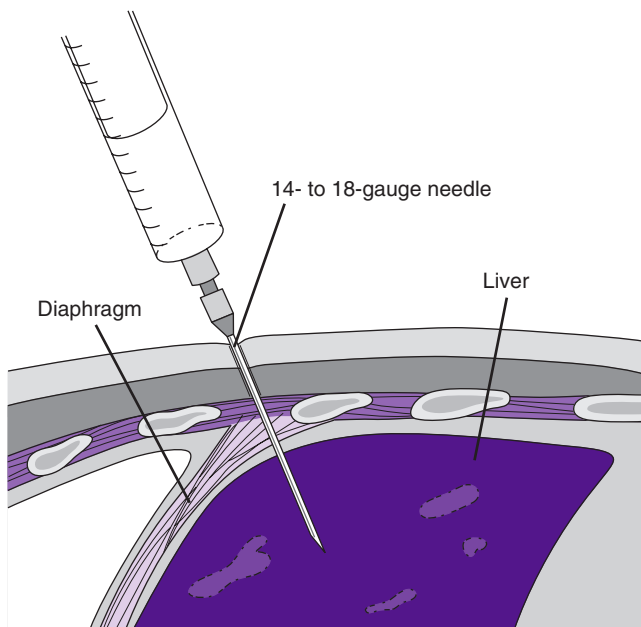

**FIG. 26** Liver biopsy. Percutaneous liver biopsy requires the patient's cooperation. The patient must be able to lie quietly and hold his or her breath after exhaling.

## 572 liver biopsy

- Patients with suspected primary or metastatic tumor
- Patients with unexplained jaundice
- Patients with suspected hepatitis
- Patients with suspected infiltrative diseases

The biopsy may be performed by a *blind stick* or directed with the use of computed tomography (CT), magnetic resonance imaging (MRI), or ultrasound.

### Contraindications

- Uncooperative patients who cannot remain still and hold their breath during sustained exhalation
- Patients with impaired hemostasis
- Patients with anemia who could not tolerate blood loss associated with puncture of an intrahepatic blood vessel
- Patients with infections in the right pleural space or right upper quadrant
- Patients with obstructive jaundice
- Patients with a hemangioma because bleeding after a biopsy may be severe

### Potential complications

- Hemorrhage caused by inadvertent puncture of a blood vessel
- Chemical peritonitis caused by inadvertent puncture of a bile duct, with subsequent leakage of bile
- Pneumothorax caused by improper placement of the biopsy needle into the adjacent chest cavity

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. Many patients are apprehensive about it.
- Obtain informed consent.
- Ensure that all coagulation tests are normal.
- PT** Instruct the patient to keep NPO after midnight on the day of the test.
- Administer any sedative medications as ordered.

#### During

- Note the following procedural steps:
  1. The patient is placed in the supine or left lateral position.
  2. The skin area used for puncture is anesthetized locally.
  3. The patient is asked to exhale and hold the exhalation. This causes the liver to descend and reduces the possibility of a pneumothorax.



4. During the patient's sustained exhalation, the physician rapidly introduces the biopsy needle into the liver and obtains liver tissue. Occasionally, the biopsy needle is inserted under CT guidance. This is especially useful when tissue from a specific area of the liver is needed.

5. The needle is withdrawn from the liver.

- Note that this test is performed by a physician in approximately 15 minutes.
- PT** Inform the patient that he or she may have minor discomfort during injection of the local anesthetic and needle insertion.

### After

- Place the tissue sample into a specimen bottle containing formalin and send it to the pathology department.
  - Apply a small dressing over the needle insertion site.
  - Place the patient on his or her right side for approximately 1 to 2 hours. In this position, the liver capsule is compressed against the chest wall, thereby decreasing the risk of hemorrhage or bile leak.
  - Assess the patient's vital signs frequently for evidence of hemorrhage and peritonitis.
  - Evaluate the rate, rhythm, and depth of respirations. Report chest pain and signs of dyspnea, cyanosis, and restlessness, which may be indicative of pneumothorax.
- PT** Tell the patient to avoid coughing or straining, which may cause increased intraabdominal pressure.

### Abnormal findings

Benign tumor

Malignant tumor (primary or metastatic)

Abscess

Cyst

Hepatitis

Infiltrative diseases (e.g., amyloidosis, hemochromatosis, cirrhosis)

notes

**liver/spleen scanning** (Liver scanning)**Type of test** Nuclear scan**Normal findings**

Normal size, shape, and position of the liver and spleen

**Test explanation and related physiology**

This radionuclide procedure is used to outline and detect structural changes of the liver and spleen. Single-photon emission computed tomography (SPECT) has significantly improved the quality and accuracy of liver scanning. With the use of radioactive carbon, nitrogen, or oxygen, anatomic and biochemical changes can be visualized within the liver. This method of liver scanning is called positron emission tomography scanning (see p. 719).

Nuclear liver scan can also identify portal hypertension. Normally, most of the radionuclide administered during a liver scan is taken up by the liver. If the liver-to-spleen ratio is reversed (i.e., the spleen takes up more of the radionuclide), reversal of hepatic blood flow exists as a result of portal hypertension.

**Contraindications**

- Patients who are pregnant or lactating unless the benefits of testing outweigh the risk of damage to the fetus or infant

**Interfering factors**

- Barium in the gastrointestinal tract overlying the liver or spleen will produce defects on the scan that may be mistaken for masses.

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Tell the patient that no fasting or premedication is required.
- PT** Assure the patient that he or she will not be exposed to large amounts of radiation.

**During**

- Note the following procedural steps:
  1. The patient is taken to nuclear medicine, where the radionuclide is administered intravenously.
  2. Thirty minutes after injection, a gamma ray detector is placed over the right upper quadrant of the abdomen.

3. The patient is placed in supine, lateral, and prone positions so that all surfaces of the liver can be visualized.
  4. The radionuclide image is recorded digitally or on an analog film.
- Note that this procedure is performed by a trained technologist in approximately 1 hour. A physician trained in nuclear medicine interprets the results.
- PT** Tell the patient that the only discomfort associated with this procedure is the (intravenous) injection of the radionuclide.

### After

- PT** Because only tracer doses of radioisotopes are used, inform the patient that no radiation precautions are needed.

### Abnormal findings

Primary or metastatic tumor of the liver or spleen

Abscess of the liver or spleen

Hematoma of the liver or spleen

Hepatic or splenic cyst

Hemangioma

Lacerations of the liver or spleen

Infiltrative processes (e.g., sarcoidosis, amyloidosis, tuberculosis, or granuloma of the liver or spleen)

Cirrhosis

Portal hypertension

Accessory spleen

Splenic infarction

### notes

**lumbar puncture and cerebrospinal fluid examination** (LP and CSF examination, Spinal tap, Cerebrospinal fluid analysis)

**Type of test** Fluid analysis

**Normal findings**

Pressure: < 20 cm H<sub>2</sub>O

Color: clear and colorless

Blood: none

Cells:

RBC count: 0

WBC count:

Total

Neonate: 0-30 cells/ $\mu$ L

1-5 years: 0-20 cells/ $\mu$ L

6-18 years: 0-10 cells/ $\mu$ L

Adult: 0-5 cells/ $\mu$ L

Differential

Neutrophils: 0%-6%

Lymphocytes: 40%-80%

Monocytes: 15%-45%

Culture and sensitivity: no organisms present

Protein: 15-45 mg/dL CSF (up to 70 mg/dL in elderly adults and children)

Protein electrophoresis

Prealbumin: 2%-7%

Albumin: 56%-76%

Alpha1 globulin: 2%-7%

Alpha2 globulin: 4%-12%

Beta globulin: 8%-18%

Gamma globulin: 3%-12%

Oligoclonal bands: none

Immunoglobulin G (IgG): 0.0-4.5 mg/dL

Glucose: 50-75 mg/dL CSF or 60%-70% of blood glucose level

Chloride: 700-750 mg/dL

LDH:  $\leq$  40 units/L (adults),  $\leq$  70 units/L (neonates)

Lactic acid: 10-25 mg/dL

Cytology: no malignant cells

Serology for syphilis: negative

Glutamine: 6-15 mg/dL

## Test explanation and related physiology

By placing a needle in the subarachnoid space of the spinal column, one can measure the pressure of that space and obtain CSF for examination and diagnosis. Lumbar puncture (LP) may also be used therapeutically to inject therapeutic or diagnostic agents and to administer spinal anesthetics. Furthermore, LP may be used to reduce intracranial pressure (ICP) in patients with normal pressure hydrocephalus or in patients with pseudotumor cerebri.

Examination of the CSF includes evaluation for the presence of blood, bacteria, and malignant cells along with quantification of the amount of glucose and protein present. Color is noted, and various other tests, such as a serologic test for syphilis (see p. 859), are performed.

### Pressure

By attaching a sterile manometer to the needle used for LP, the pressure within the subarachnoid space can be measured. A pressure greater than 20 cm H<sub>2</sub>O is considered abnormal and indicative of increased spinal pressure. Because the subarachnoid space surrounding the brain is freely connected to the subarachnoid space of the spinal cord, any increase in ICP will be directly reflected as an increase at the lumbar site. Tumors, infection, hydrocephalus, and intracranial bleeding can cause increased intracranial and spinal pressure. ICP is related to the volume of CSF fluid. Also, because the cranial venous sinuses are connected to the jugular veins, obstruction of those veins or the superior vena cava increases ICP.

### Color

Normal CSF is clear and colorless. Color differences can occur with hyperbilirubinemia, hypercarotenemia, melanoma, or elevated proteins. A cloudy appearance may indicate an increase in the WBC count or protein. A red tinge to the CSF indicates the presence of blood.

### Blood

Normally, CSF contains no blood. Blood may be present because of bleeding into the subarachnoid space or because the needle used in the LP has inadvertently penetrated a blood vessel on the way into the subarachnoid space. With a *traumatic puncture*, the blood within the CSF will clot. No clotting occurs in a patient with subarachnoid hemorrhage. Also, with a traumatic tap, the fluid clears toward the end of the procedure as successive CSF samples are obtained. This clearing does not occur with a subarachnoid hemorrhage.

## 578 lumbar puncture and cerebrospinal fluid examination

### Cells

The number of red blood cells (RBCs) is merely an indication of the amount of blood present within the CSF. Except for a few lymphocytes, the presence of WBCs in the CSF is abnormal. The presence of polymorphonuclear leukocytes (neutrophils) is indicative of bacterial meningitis or cerebral abscess. When mononuclear leukocytes are present, viral or tubercular meningitis or encephalitis is suspected. Leukemia or other primary or metastatic malignant tumors may cause elevated WBCs.

### Culture and sensitivity

The organisms that cause meningitis or brain abscess can be cultured from the CSF. Organisms found also may include atypical bacteria, fungi, or *Mycobacterium tuberculosis*. A Gram stain (see p. 841) of the CSF may give the clinician preliminary information about the causative infectious agent.

### Protein

Normally, very little protein is found in CSF because proteins are large molecules that do not cross the blood-brain barrier. Normally, the proportion of albumin to globulin is higher in CSF than in blood plasma (see p. 746) because albumin is smaller than globulin and can pass more easily through the blood-brain barrier. Such diseases as meningitis, encephalitis, and myelitis can alter the permeability of the blood-brain barrier, allowing protein to leak into the CSF. Furthermore, central nervous system (CNS) tumors may produce and secrete protein into the CSF.

CSF protein electrophoresis is very important in the diagnosis of CNS diseases. Patients with multiple sclerosis (MS), neurosyphilis, or other immunogenic degenerative central neurologic diseases have elevated immunoglobulins in their CSF. The detection of *oligoclonal gamma globulin bands* is highly suggestive of inflammatory and autoimmune diseases of the CNS, especially MS.

### Glucose

The glucose level is decreased when bacteria, inflammatory cells, or tumor cells are present. A blood sample for glucose (see p. 462) is usually drawn before the spinal tap is performed. A CSF glucose level of less than 60% of the blood glucose may indicate meningitis or neoplasm.

### Chloride

The chloride concentration in CSF may be decreased in patients with meningeal infections, tubercular meningitis, and conditions of low blood chloride levels.

**Lactic dehydrogenase**

Quantification of LDH (specifically fractions 4 and 5; see p. 553) is helpful in diagnosing bacterial meningitis. The sources of LDH are the neutrophils that fight the invading bacteria. When the LDH level is elevated, infection or inflammation is suspected. The elevated WBC count associated with CNS leukemia is also associated with elevated LDH levels. The nerve tissue in the CNS is also high in LDH (isoenzymes 1 and 2). Therefore disease directly affecting the brain or spinal cord (e.g., stroke) are associated with elevated LDH levels.

**Lactic acid**

Elevated levels indicate anaerobic metabolism associated with decreased oxygenation of the brain. CSF lactic acid is increased in both bacterial and fungal meningitis but not in viral meningitis.

**Cytology**

Examination of cells found in the CSF can determine whether they are malignant. Tumors in the CNS may shed cells from their surface. These cells can float freely in CSF.

**Tumor markers**

Increased levels of tumor markers (see p. 194) (e.g., carcinoembryonic antigen, alpha-fetoprotein, and human chorionic gonadotropin) may indicate metastatic tumor.

**Serology for syphilis**

Latent syphilis is diagnosed by performing one of many presently available serologic tests on CSF. These include the Wasserman test, the Venereal Disease Research Laboratory test, and the fluorescent treponemal antibody (FTA) test (see p. 859). When test results are positive, the diagnosis of neurosyphilis is made, and appropriate antibiotic therapy is initiated.

**Glutamine**

Elevated glutamine levels are helpful in the detection and evaluation of hepatic encephalopathy and hepatic coma. The glutamine is made by increased levels of ammonia, which are commonly associated with liver failure (see serum ammonia, p. 47).

**C-reactive protein**

As noted on p. 295, C-reactive protein (CRP) is a nonspecific, acute-phase reactant used in the diagnosis of bacterial infections and inflammatory disorders. Elevated CSF levels of CRP have been useful in the diagnosis of bacterial meningitis. Failure to find elevated CSF levels of CRP appears to be strong evidence against bacterial meningitis.

## 580 lumbar puncture and cerebrospinal fluid examination

### Contraindications

- Patients with increased ICP. The LP may induce cerebral or cerebellar herniation through the foramen magnum.
- Patients who are anticoagulated. Because of the risk of epidural hematoma, these patients need to have normal coagulation function before this procedure.
- Patients who have severe degenerative vertebral joint disease. It is very difficult to pass the needle through the degenerated arthritic interspinal space.
- Patients with infection near the LP site. Meningitis can result from contamination of CSF.

### Potential complications

- Persistent CSF leak, causing severe headache
- Puncture of subcutaneous blood vessel during the procedure
- Introduction of bacteria causing meningitis
- Herniation of the brain through the tentorium cerebelli or herniation of the cerebellum through the foramen magnum
- Inadvertent puncture of the spinal cord caused by inappropriately high puncture of the spinal canal
- Puncture of the aorta or vena cava, causing serious retroperitoneal hemorrhage
- Transient back pain and pain or paresthesia in the legs
- Transient postural headache (worse when standing)

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. Many patients have misconceptions regarding LP.
- Obtain informed consent if required by the institution.
- Perform a baseline neurologic assessment of the legs.
- PT** Tell the patient that no fasting or sedation is required.
- PT** Instruct the patient to empty the bladder and bowels.
- PT** Explain to the patient that he or she must lie very still throughout this procedure. Movement may cause injury.

#### During

- Note the following procedural steps:
  1. This study can be easily performed at the bedside. The patient is usually placed in the lateral decubitus (fetal) position ([Figure 27](#)).
  2. The patient is instructed to clasp the hands on the knees to maintain this position. (A sitting position also may be used.)



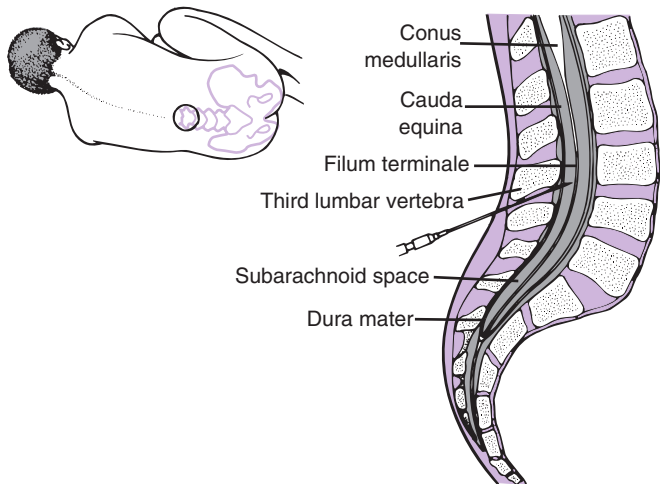


FIG. 27 Patient position for lumbar puncture.

3. A local anesthetic is injected into the skin and subcutaneous tissues after the site has been aseptically cleaned.
  4. A spinal needle containing an inner obturator is placed through the skin and into the spinal canal.
  5. The subarachnoid space is entered.
  6. The insert (obturator) is removed, and CSF can be seen slowly dripping from the needle.
  7. The needle is attached to a sterile manometer, and the pressure (opening pressure) is recorded.
  8. Before the pressure reading is taken, the patient is asked to relax and straighten the legs to reduce the intraabdominal pressure.
  9. Three sterile test tubes are filled with 5 to 10 mL of CSF.
  10. The pressure (closing pressure) is measured.
- Note that if blockage in CSF circulation in the spinal subarachnoid space is suspected, a *Queckenstedt-Stookey* test may be performed. For this test, the jugular vein is occluded either manually by digital pressure or by a medium-sized blood pressure cuff inflated to approximately 20 mm Hg. Within 10 seconds after jugular occlusion, CSF pressure should increase from 15 to 40 cm H<sub>2</sub>O and then promptly return to normal within 10 seconds after release of the pressure. A sluggish rise

## 582 lumbar puncture and cerebrospinal fluid examination

or fall of CSF pressure suggests partial blockage of CSF circulation. No rise after 10 seconds suggests a complete obstruction within the spinal canal.

- Note that this procedure is performed by a physician in approximately 20 minutes.
- PT** Inform the patient that this procedure is described as uncomfortable or painful by most patients.

### After

- Apply digital pressure and a dressing to the puncture site.
  - Place the patient in the prone position with a pillow under the abdomen to increase the intraabdominal pressure, which will indirectly increase the pressure in the tissues surrounding the spinal cord. This retards continued CSF flow from the spinal canal.
- PT** Encourage the patient to drink increased amounts of fluid with a straw to replace the CSF removed during LP.
- PT** Usually keep the patient in a reclining position for 1 hour or up to several hours to avoid the discomfort of potential post-puncture spinal headache. Instruct the patient to turn from side to side as long as the head is not raised.
- Label and number the specimen jars appropriately and deliver them immediately to the laboratory after the test.
- PT** Instruct the patient to report numbness, tingling, and movement of the extremities; pain at the injection site; drainage of blood or CSF at the injection site; and the inability to void.

### Abnormal findings

Brain neoplasm  
Spinal cord neoplasm  
Cerebral hemorrhage  
Encephalitis  
Myelitis  
Tumor  
Neurosyphilis  
Autoimmune disorder  
Hepatic encephalopathy  
Coma  
Meningitis  
Viral or tubercular meningitis  
Cerebral abscess  
Degenerative cord or brain disease

Multiple sclerosis  
Acute demyelinating polyneuropathy  
Subarachnoid bleeding  
Reye syndrome  
Metastatic tumor

**notes**

## lung biopsy

**Type of test** Microscopic examination of tissue

### Normal findings

No evidence of pathology

### Test explanation and related physiology

This invasive procedure is used to obtain a specimen of pulmonary tissue for a histologic examination by using either an open or a closed technique. The *open method* involves a limited thoracotomy. The *closed technique* includes methods such as transbronchial lung biopsy, transbronchial needle aspiration biopsy, transcatheter bronchial brushing, percutaneous needle biopsy, and video-assisted thoracotomy (see p. 870).

Lung biopsy is indicated to determine the pathology of pulmonary parenchymal disease. Carcinomas, granulomas, infections, and sarcoidosis can be diagnosed with this procedure. The procedure is also useful in detecting environmental exposures, infections, and familial disease.

### Contraindications

- Patients with bullae or cysts of the lung
- Patients with suspected vascular anomalies of the lung
- Patients with bleeding abnormalities
- Patients with pulmonary hypertension
- Patients with respiratory insufficiency

### Potential complications

- Pneumothorax
- Pulmonary hemorrhage
- Empyema

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Ensure that informed and signed consent is obtained.

**PT** Instruct the patient that fasting is usually ordered. The patient may be kept NPO (nothing by mouth) after midnight on the day of the test.

**PT** Instruct the patient to remain still during the lung biopsy.

#### During


- Note that the patient's position depends on the method used and that the histologic lung specimen may be obtained by several different methods.

***Transbronchial lung biopsy***

- This technique is performed via flexible fiberoptic bronchoscopy by using cutting forceps.
- Fluoroscopy is used to ensure proper opening and positioning of the forceps on the lesions.
- Fluoroscopy also permits visualization of the tug of the lung as the specimen is removed.

***Transbronchial needle aspiration***

- The needle is inserted through the bronchoscope and into the tumor or desired area, where aspiration is performed with the attached syringe (Figure 28).
- The needle is retracted within its sheath, and the entire catheter is withdrawn from the fiberoptic scope.

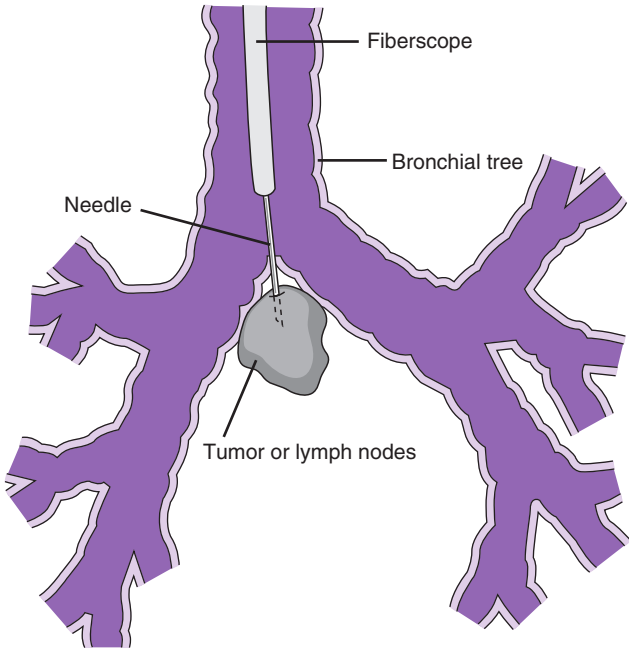

**FIG. 28** Transbronchial needle biopsy. The diagram shows a transbronchial needle penetrating the bronchial wall and entering a mass of subcarinal lymph nodes or tumor.

*Transbronchial brushing*

- A small brush is moved back and forth over the suspicious area in the bronchus or its branches.
- The cells adhere to the brush, which is then removed and used to make microscopic slides.

*Percutaneous needle biopsy*

- In this method of obtaining a closed specimen, the biopsy is obtained after using fluoroscopic radiography or CT scan.
- The procedure is carried out by using a cutting needle or by aspiration with a spinal-type needle to obtain a specimen.

*Open lung biopsy*

- The patient is taken to the operating room, and general anesthesia is provided.
- An incision is made into the chest wall.
- After a piece of lung tissue is removed, the lung is sutured.
- Chest tube drainage is used for approximately 24 hours.

*Thoracoscopic lung biopsy*

- The lung is collapsed with a double-lumen endotracheal tube placed during induction of general anesthesia.
  - With the use of a thoracoscope, the lung is grasped, and a piece is cut off with the use of a cutting or stapling device.
  - The scope and trocars are removed, and a small chest tube is left in place. The tiny incisions are closed.
  - Note that this procedure is performed by a surgeon, radiologist, or pulmonologist in 30 to 60 minutes.
  - During the lung biopsy procedure, assess the patient carefully for signs of respiratory distress.
- PT** Tell the patient that most patients describe this procedure as painful.

**After**

- Place biopsy specimens in appropriate containers for histologic and microbiologic examination.
- Observe the patient's vital signs frequently for signs of bleeding and shortness of breath.
- Assess the patient's breath sounds and report any decrease.
- Obtain a chest x-ray image to check for complications.
- Observe the patient for signs of pneumothorax (e.g., dyspnea, tachypnea, decrease in breath sounds, anxiety, restlessness).

## Abnormal findings

Carcinoma

Granuloma

Exposure lung diseases (e.g., black lung, asbestosis)

Sarcoidosis

Infection

notes

**lung cancer molecular testing** (AKT1, ALK, BRAF, Epidermal growth factor [EGFR], HER2, KRAS, MEK1, MET, NRAS, PIK3CA, RET, and ROS1)

**Type of test** Microscopic examination

### Normal findings

Absent genetic mutations

### Test explanation and related physiology

Treatment for lung cancer has been empiric and based on histology of the tumor (i.e., adenocarcinoma, squamous, small cell). Subsets of non–small cell lung cancers (NSCLC) can be further defined at the molecular level by the presence of mutations that may occur in oncogenes (e.g., AKT1, ALK, BRAF, EGFR, HER2, KRAS, MEK1, MET, NRAS, PIK3CA, RET, and ROS1). Mutations in these oncogenes lead to activation of signaling proteins that induce and sustain tumor growth and therapy resistance. These mutations are not commonly found concurrently in the same tumor. Mutations can be found in all NSCLC histologies (including adenocarcinoma, squamous cell carcinoma [SCC], and large cell carcinoma) and in smokers and nonsmokers, alike.

Importantly, targeted drugs are currently available and still more are being developed that can inhibit the tumor development of these proteins and thereby inhibit tumor growth. Most significantly, the therapeutic importance of genes that encode pharmacologically targetable tyrosine kinases involved in growth factor receptor signaling, epidermal growth factor receptor (EGFR), and anaplastic lymphoma kinase (ALK) has changed the way these cancers are diagnosed and treated. At present, expert panels of oncologists and laboratory pathologists are attempting to develop consensus on the proper use of these tests in the treatment of local and metastatic lung cancer. Their use must take into account their cost of testing.

### Procedure and patient care

#### Before

- PT** Explain the significance of the diagnostic data available.
- PT** Explain the benefits of molecular testing in helping the physician and the patient make appropriate decisions regarding the use of surgery.



- Provide the patient with emotional support through the diagnostic period.
- Ensure that the patient's insurance will cover this expensive testing.

### During

- The pathologists will send a fresh or paraffin-embedded tumor specimen to a centralized laboratory.
- Results will be available in about 2 weeks.

### After

- PT** Provide education and support to patients as they evaluate their results.

### Abnormal findings

Known mutation of respective gene.

notes

**lung scan** (Ventilation/perfusion scanning [VPS], Pulmonary scintiphotography, V/Q scan)**Type of test** Nuclear scan**Normal findings**

Diffuse and homogeneous uptake of nuclear material by the lungs

**Test explanation and related physiology**

This nuclear medicine procedure is used to identify defects in blood *perfusion* of the lung in patients with suspected pulmonary embolism. A homogeneous uptake of particles that fills the entire pulmonary vasculature conclusively rules out pulmonary embolism. If a defect in an otherwise smooth and diffusely homogeneous pattern is seen, a perfusion abnormality exists. This can indicate pulmonary embolism. Unfortunately, many other serious pulmonary parenchymal lesions (e.g., pneumonia, pleural fluid, emphysematous bullae) also cause a defect in pulmonary blood perfusion. Therefore although the scan may be sensitive, it is not specific because many different pathologic conditions can cause the same abnormal results.

The chest x-ray image aids in assessing the perfusion scan because a defect on the perfusion scan seen in the same area as an abnormality on the chest x-ray image does not indicate pulmonary embolism. Rather, the defect may represent pneumonia, atelectasis, effusion, and so on. However, when a perfusion defect occurs in an area of the lung that is normal on a chest x-ray study, pulmonary embolus is likely.

Specificity of a perfusion scan also can be enhanced by performance of a ventilation scan, which detects parenchymal abnormalities in ventilation (e.g., pneumonia, pleural fluid, emphysematous bullae). The *ventilation scan* reflects the patency of the pulmonary airways by using xenon gas or Tc-diethylenetriamine pentaacetic acid (DTPA) as an aerosol. When vascular obstruction (embolism) is present by perfusion scan, ventilation scans demonstrate a normal wash-in and wash-out of radioactivity from the embolized lung area. If parenchymal disease (e.g., pneumonia) is responsible for the perfusion abnormality, however, wash-in or wash-out will be abnormal. Therefore whereas the mismatch of perfusion and ventilation findings is characteristic of embolic disorders, the match is indicative of parenchymal disease. When ventilation and perfusion scans are performed synchronously, this is called a *ventilation/perfusion ( $\dot{V}/\dot{Q}$ ) scan*. Most nuclear

physicians read the lung scan as one of several categories: negative for PE, low probability of PE, high probability of PE, or positive for PE.

### Contraindications

- Patients who are pregnant unless the benefits outweigh the risks

### Interfering factors

**PT** Pulmonary parenchymal problems (e.g., pneumonia, emphysema, pleural effusion, tumors) will give the picture of a perfusion defect and simulate pulmonary embolism.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.

- Obtain informed consent if required by the institution.

**PT** Assure the patient that he or she will not be exposed to large amounts of radioactivity because only tracer doses of isotopes are used.

**PT** Tell the patient that no fasting is required.

- Note that a recent chest x-ray image should be available.

**PT** Instruct the patient to remove jewelry around the chest.

#### During

- Note the following procedural steps for a V/Q scan:
  1. The unsedated, nonfasting patient is taken to the nuclear medicine department.

#### *Ventilation scan*

2. The patient breathes the tracer through a face mask with a mouthpiece.
3. Ventilation scans require a cooperative patient.
4. Technetium-diethylenetriaminepentaacetic acid images are usually done before perfusion images and require patient cooperation with deep breathing and appropriate use of breathing equipment to prevent contamination.

#### *Perfusion scan*

5. The patient is given a peripheral IV injection of radionuclide-tagged macroaggregated albumin (MAA).
6. While the patient lies in the appropriate position, a gamma ray detector is passed over the patient and records radionuclide uptake on film.
7. The patient is placed in the supine, prone, and various lateral positions, which allow for anterior, posterior, lateral, and oblique views, respectively.

## 592 lung scan

8. The results are interpreted by a physician trained in diagnostic nuclear medicine.
- Note that this test is usually performed by a technologist in approximately 30 minutes.
- PT** Tell the patient that no discomfort is associated with this test other than the peripheral venipuncture.

### After

- Apply pressure to the venipuncture site.
- PT** Inform the patient that no radiation precautions are necessary.

### Abnormal findings

Pulmonary embolism

Pneumonia

Tuberculosis

Emphysema

Tumor

Asthma

Atelectasis

Bronchitis

Chronic obstructive pulmonary disease

### notes

## luteinizing hormone (LH assay, Lutropin) and follicle-stimulating hormone (FSH) assay

**Type of test** Blood

### Normal findings

	LH (IU/L)	FSH (IU/L)
<i>Adult</i>		
Male	1.24-7.8	1.42-15.4
Female		
Follicular phase	1.68-15	1.37-9.9
Ovulatory peak	21.9-56.6	6.17-17.2
Luteal phase	0.61-16.3	1.09-9.2
Postmenopause	14.2-52.3	19.3-100.6
<i>Child</i>		
(age 1-10 years)		
Male	0.04-3.6	0.3-4.6
Female	0.03-3.9	0.68-6.7

(Values may vary, depending on laboratory method.)

### Test explanation and related physiology

LH and FSH are glycoproteins that are produced in the anterior pituitary gland. These two hormones then act on the ovary or testes. Spot urine tests for LH have become very useful in the evaluation and treatment of infertility. Because LH is rapidly excreted into the urine, the plasma LH surge that precedes ovulation by 24 hours can be recognized quickly and easily. This is used to indicate the period when a woman is the most fertile. The best time to obtain a urine specimen is between 11 AM and 3 PM. Usually, the woman begins to test her urine on the 10th day after the onset of her menses and continues to do so daily. Home kits in which a color change as an endpoint are used to make this process more convenient.

These hormones are used in the evaluation of infertility. Performing an LH assay is an easy way to determine whether ovulation has occurred. An LH surge in blood levels indicates that ovulation has taken place. Daily samples of serum LH around the woman's midcycle can detect the LH surge, which is believed to occur on the day of maximal fertility.

A more accurate method to evaluate ovarian function is by ovarian reserve testing. A commonly used test is the *Clomiphene (Clomid) Challenge Test (CCCT)*. In this test, clomiphene is

## 594 luteinizing hormone

administered on days 5 to 9 of a woman's menstrual cycle. FSH levels are expected to rise significantly and then fall considerably by day 10 of the menstrual cycle. Persistently elevated levels indicate poor ovarian reserve and predict a reduced chance of pregnancy.

### Interfering factors

- Human chorionic gonadotropin (hCG) and thyroid-stimulating hormone may interfere with some immunoassay methods. Patients with hCG-producing tumors and hypothyroid patients may have falsely high LH levels.
- ▼ Drugs that may *increase* LH levels include anticonvulsants, clomiphene, naloxone, and spironolactone.
- ▼ Drugs that may *decrease* LH levels include digoxin, estrogens, oral contraceptives, phenothiazines, progesterones, and testosterone.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Note that the patient may also perform LH assays at home using a home urine test or a 24-hour urine test.
- Indicate the date of the last menstrual period on the laboratory slip. Note if the woman is postmenopausal.

### Abnormal findings

#### ▲ Increased levels

Menopause  
Ovarian dysgenesis (Turner syndrome)  
Testicular dysgenesis (Klinefelter syndrome)  
Castration  
Anorchia  
Hypogonadism  
Polycystic ovaries  
Complete testicular feminization syndrome  
Precocious puberty  
Pituitary adenoma

#### ▼ Decreased levels

Pituitary failure  
Hypothalamic failure  
Stress  
Anorexia nervosa  
Malnutrition

## lyme disease test

**Type of test** Blood

### Normal findings

*Borrelia burgdorferi* antibody Enzyme Immunoassay (Lyme index value)

< 0.9 = negative

0.91-1.09 = equivocal

> 1.1 = positive

### Western blot

≥ 5 different IgG antibodies reactive = positive

≥ 2 different IgM antibodies reactive = positive

PCR: negative

CSF: negative

### Test explanation and related physiology

Lyme disease is caused by a spirochete called *B. burgdorferi*. This is the most common tickborne disease. The spirochete is spread by a bite from a black-legged tick (*Ixodes pacificus*) or deer tick (*Ixodes scapularis*).

Cultures of the distinctive primary skin lesion (erythema chronicum migrans) can isolate the spirochete in half of the cases. However, it is difficult to culture and takes a long time to grow. Cultures of the blood or CSF are even less helpful. Currently, screening serologic studies are performed for the detection of Lyme disease. ELISA is the best diagnostic test for Lyme disease. This test determines titers of specific IgM and specific IgG antibodies to the *B. burgdorferi* spirochete. Levels of specific IgM antibody peak during the third to sixth week after disease onset and then gradually decline. Titers of specific IgG antibodies are generally low during the first several weeks of illness, reach maximal levels in 4 to 6 months, and often remain elevated for years.

Lyme disease can be confused with various viral infections. In these patients, a single titer of specific IgM antibody may suggest the correct diagnosis. Acute and convalescent sera can be tested to verify the diagnosis with a significant rise in positive antibody titers. The Food and Drug Administration recommends that all samples with positive or equivocal results in the *B. burgdorferi* antibody ELISA (screening) should be tested by Western blot. Positive or equivocal ELISA screening test results should not be interpreted as truly positive until verified with a confirmatory Western blot assay. The Western blot antibody assay can

## 596 lyme disease test

identify specifically the IgG or the IgM antibody. The Western blot assay is considered positive for IgG if 5 or more of the 10 significant electrophoretic bands are considered positive for *B. burgdorferi*-specific IgG antibody. The Western blot IgM antibody assay is considered positive if two or more of three significant electrophoretic bands are considered positive for *B. burgdorferi* IgM antibody. However, the screening test or Western blot for *B. burgdorferi* antibodies may be falsely negative in early stages of Lyme disease, including the period when erythema migrans is apparent.

### Interfering factors

- Previous infection with *B. burgdorferi* can cause positive serologic testing results. These patients may no longer have Lyme disease.
- Other spirochete diseases (syphilis or leptospirosis) can cause false-positive results.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

Lyme disease

notes



## magnesium

**Type of test** Blood

### Normal findings

Adult: 1.3-2.1 mEq/L or 0.65-1.05 mmol/L (SI units)

Child: 1.4-1.7 mEq/L

Newborn: 1.4-2 mEq/L

### Possible critical values

< 0.5 mEq/L or > 3 mEq/L

### Test explanation and related physiology

Most organ functions also depend on magnesium. It is especially important to monitor magnesium levels in cardiac patients. Low magnesium may increase cardiac irritability and aggravate cardiac arrhythmias. Hypermagnesemia retards neuromuscular conduction and is demonstrated as cardiac conduction slowing (widened PR and Q-T intervals, with wide QRS).

As intracellular elements, potassium, magnesium, and calcium (in order of quantity) are intimately tied together in maintaining a neutral intracellular electrical charge. That is why it is hard to maintain a normal potassium level when a patient has low magnesium blood levels.

Magnesium deficiency occurs in patients who are malnourished. Toxemia of pregnancy is also believed to be associated with reduced magnesium levels. Magnesium testing is used to monitor magnesium replacement therapy.

Increased magnesium levels most commonly are associated with ingestion of magnesium-containing antacids. Because most of the serum magnesium is excreted by the kidneys, chronic renal diseases cause elevated magnesium levels.

### Interfering factors

- Hemolysis should be avoided when collecting this specimen.
- ✚ Drugs that *increase* magnesium levels include aminoglycoside antibiotics, antacids, calcium-containing medications, laxatives, lithium, loop diuretics, and thyroid medication.
- ✚ Drugs that *decrease* magnesium levels include some antibiotics, diuretics, and insulin.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

▲ Increased levels

- Renal insufficiency
- Uncontrolled diabetes
- Addison disease
- Hypothyroidism
- Ingestion of magnesium-containing antacids or salts

▼ Decreased levels

- Malnutrition
- Malabsorption
- Hypoparathyroidism
- Alcoholism
- Chronic renal disease
- Diabetic acidosis

notes

## magnetic resonance imaging (MRI, Nuclear magnetic resonance imaging [NMRI])

**Type of test** Magnetic field study

### Normal findings

No evidence of pathology

### Test explanation and related physiology

MRI is a noninvasive diagnostic scanning technique that provides valuable information about the body's anatomy by placing the patient in a magnetic field. MRI does not require exposure to ionizing radiation.

MRI is useful in the evaluation of the following areas:

- Head and surrounding structures
- Spinal cord and surrounding structures
- Face and surrounding structures
- Neck
- Mediastinum
- Heart and great vessels
- Liver and biliary tree
- Kidney
- Prostate
- Bones and joints
- Breast
- Extremities and soft tissues
- Pancreas

*MRI of the brain and meninges* is particularly accurate in identifying benign and malignant neoplasms. It is able to identify and quantify brain edema, ventricular compression, hydrocephalus, and brain herniation. Intracranial hemorrhage can also be seen on MRI. *Magnetic resonance spectroscopy (MRS)* is a noninvasive procedure that generates high-resolution clinical images based on the distribution of chemicals in the body. This is particularly useful in the brain, in which certain chemical metabolites will enhance the image of a high-grade malignancy. MRS has also been used to assess chemical abnormalities in the brain associated with HIV infection without having to perform a brain biopsy. This procedure has been used in a wide variety of disorders, including stroke, head injury, coma, Alzheimer disease, and multiple sclerosis.

MRI has revolutionized the practice of orthopedic surgery. It is particularly helpful in the determination of anatomic changes in muscle and joints (particularly knee and shoulder).

## 600 magnetic resonance imaging

*Magnetic resonance angiography (MRA)* is a noninvasive procedure for viewing possible blockages in arteries. MRA has been useful in evaluation of the cervical carotid artery and large-caliber intracranial arterial and venous structures. Cardiac abnormalities, aortic aneurysm, and anatomic variants can be identified. This procedure also has proved useful in the noninvasive detection of intracranial aneurysms and vascular malformations, especially in renal artery stenosis.

*MRI of the breast* is more sensitive than mammography or ultrasonography of the breast. MRI of the breast is used for accurate localized staging of breast cancer by demonstrating an excellent three-dimensional image of a cancer and high sensitivity for other smaller synchronously occurring breast cancers that may be missed on mammography. MRI of the breast is helpful for preoperative surgical staging and for the identification of postoperative positive margins. This study is particularly helpful in differentiating postoperative scar tissue from breast cancer recurrence. MRI of the breast is the most accurate method of determining fracture of a breast implant. With the addition of a needle-guiding system to the MRI, breast tumors can be nonoperatively and accurately localized and biopsied.

Significant improvement in *MRI of the heart* and great vessels has moved this noninvasive diagnostic procedure into the mainstream of clinical cardiology. Cardiac MRI already is considered the procedure of choice in the evaluation of pericardial disease and intracardiac and pericardiac masses, for imaging the right ventricle and pulmonary vessels, and for assessing many forms of congenital heart disease, especially after corrective surgery. The ventricle size, shape, and blood volumes can be evaluated. Cardiac valvular abnormalities, cardiac septal defects, and suspected intracardiac or pericardiac masses or thrombi can be identified. Pericardial disease (e.g., pericarditis or effusion) is easily identified. Ventricular muscle changes from ischemia or infarction can be determined. Finally, advanced MRI techniques are able to evaluate the coronary vessels directly.

*Phase-contrast magnetic resonance imaging (PC-MRI)* of the heart quantifies velocity and blood flow in the great arteries. Measurements of blood flow in the aorta and pulmonary trunk produce a wealth of information, including cardiac outputs of the left and right ventricles, regurgitant volumes and fraction of the aortic and pulmonary valves, and shunt ratio. Stress cardiac MRI can be performed using nitrates, dobutamine, and adenosine. When beta-blockers are added to electrocardiographic gating, cardiac volumes and images can be better portrayed.

*Magnetic resonance cholangiopancreatography (MRCP)* allows noninvasive imaging of the biliary tree, gallbladder, pancreas, and pancreatic duct. It is used to:

- Identify tumors, stones, inflammation, or infection.
- Evaluate the cause of pancreatitis.
- Help in the diagnosis of unexplained abdominal pain.
- Provide a less invasive alternative to endoscopic retrograde cholangiopancreatography (ERCP).

Indications for the use of MRCP include unsuccessful or contraindicated ERCP, patient preference for noninvasive imaging, patients considered to be at low risk of having pancreatic or biliary disease, patients in which the need for therapeutic ERCP is considered unlikely, and those with a suspected neoplastic cause for pancreatic or biliary obstruction.

*Magnetic resonance enterography (MRE)* is used to identify inflammatory bowel disease. It is also helpful in determining extraluminal bowel pathology. MRI is an effective tool in liver imaging and in the staging of known prostate cancers.

One of the most common uses is *MRI of the cervical or lumbar spine*. The main purpose of this test is to determine the cause of neck or back pain, respectively. The MRI is the most accurate test to identify herniated disc disease. Using different MRI protocols, an *MRI myelogram* can be performed when the spinal fluid appears white, and the solid tissue (discs or nerves) appears dark. Herniated discs are easily seen and graded as to their compression on the nerves. Furthermore, MRI of the spine is able to identify subtle changes associated with early infiltrating diseases such as metastatic cancer. An *upright MRI* can scan patients in any position. The upright MRI can scan patients in their positions of symptoms (e.g., pain or numbness) including weight-bearing positions such as sitting, standing, or bending. The upright MRI can provide diagnostic images of the cervical spine, the lumbar spine, and the joints over their full range of motion (e.g., cervical flexion/extension). The front-open and top-open design of the upright MRI system nearly eliminates possible claustrophobia and accommodates larger patients.

*MRI of the liver* has improved significantly with the use of a gadolinium-like contrast agent called gadoxetate (Eovist). Imaging with this agent provides extremely sharp imaging that can identify liver and biliary tumors smaller than 1 cm.

*Magnetic resonance venography (MRV)* uses magnetic resonance technology to visualize the veins. It provides images for the diagnosis of venous abnormalities. This mode of imaging is often an alternative approach for patients with contrast dye

## 602 magnetic resonance imaging

allergy. It can be used to diagnose deep vein thrombosis and is a good alternative to the more invasive contrast venography (see p. 963) and ultrasound venous duplex scan (see p. 960). MRV is also used in the diagnosis of cerebral venous thrombosis by demonstrating absence of flow in the cerebral venous sinuses. It can also help in the diagnosis of idiopathic intracranial hypertension and superior vena cava syndrome.

### Potential complications

- Gadolinium-based contrast agents (gadopentetate dimeglumine [Magnevist], gadobenate dimeglumine [MultiHance], gadodiamide [Omniscan], gadoversetamide [OptiMARK], and gadoteridol [ProHance]) have been linked to the development of nephrogenic systemic fibrosis (NSF) or nephrogenic fibrosing dermopathy (NFD). A creatinine, blood urea nitrogen, and/or estimated glomerular filtration rate (see p. 303) may be obtained, especially on adults older than the age of 60 years.

### Contraindications

- Patients who are extremely obese (usually > 300 lb)
- Patients who are confused or agitated
- Patients who are claustrophobic, if an enclosed scanner is used. This can be overcome with the administration of anti-anxiety medication.
- Patients who are unstable and require continuous life-support equipment because most monitoring equipment cannot be used inside the scanner room. Magnet-adaptive equipment is becoming available for use in MRI scanner rooms.
- Patients with implantable metal objects (e.g., pacemakers, infusion pumps, aneurysm clips, inner ear implants, and metal fragments in one or both eyes) because the magnet may move the object in the body and injure the patient. Piercings, braces, and retainers need to be removed.

### Interfering factors

- Movement during the scan may cause artifacts on MRI.
- Permanent retainers will cause an artifact on the scan.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Inform the patient that there is no exposure to radiation.
  - Obtain informed consent if required by the institution.
- PT** Tell the patient that he or she can drive without assistance after the procedure unless anti-anxiety medications are administered to treat claustrophobia.

- PT** Tell parents of young patients that they may read or talk to a child in the scanning room during the procedure. There is no risk of radiation from the procedure.
  - Assess the patient for any contraindications for testing (e.g., aneurysm clips).
- PT** If available, show the patient a picture of the scanning machine and encourage verbalization of anxieties. Some patients may experience claustrophobia. If possible, an open MRI system can be used for these patients.
- PT** Instruct the patient to remove all metal objects (e.g., dental bridges, jewelry, hair clips, belts) because they will create artifacts on the scan. The magnetic field can damage watches and credit cards. Also, movement of metal objects within the magnetic field can be detrimental to patients or staff within the field.
- PT** Tell the patient wearing a nicotine patch (or any other patch with a metallic foil backing) to remove it. These patches can become intensely hot during the MRI and cause burns.
- PT** Inform the patient that he or she will be required to remain motionless during this study. Any movement can cause artifacts on the scan.
- PT** Tell the patient that during the procedure he or she may hear a thumping sound. Earplugs are available if the patient wishes to use them.
- PT** Inform the patient that fluid or food restrictions may be required before abdominal MRI.
- PT** For comfort, instruct the patient to empty the bladder before the test.
  - For some joint evaluations, patients have their joints injected in the x-ray department via fluoroscopy with a contrast dye before their *MRI arthrogram*.

### During

- Note the following:
  1. The patient lies on a platform that slides into a tube containing the cylinder-shaped tubular magnet. However, patients can be scanned in sitting, standing, or bending positions with an upright MRI.
  2. For cardiac MRI, electrocardiography leads are applied (see p. 342).
  3. The patient is instructed to remain very still during the procedure. The patient may be asked to stop breathing for short periods of time.
  4. During the scan, the patient can talk to and hear the staff via microphone or earphones placed in the scanner.

## 604 magnetic resonance imaging

5. A contrast medium called gadolinium is a paramagnetic enhancement agent that crosses the blood–brain barrier. It is especially useful for distinguishing hypermetabolic abnormalities such as tumors. If this is to be administered, approximately 10 to 15 mL is injected in the vein. Imaging can begin shortly after the injection. No dietary restrictions are necessary before using this agent.
- Note that a qualified radiologic technologist performs this procedure in approximately 30 to 90 minutes.
- PT** Tell the patient that the only discomfort associated with this procedure may be lying still on a hard surface and a possible tingling sensation in teeth containing metal fillings. Also, an injection may be needed for administration of the contrast medium.

### After

- PT** Inform the patient that no special postprocedural care is needed.

## Abnormal findings

### Brain

Cerebral tumor  
Cerebrovascular accident  
Aneurysm  
Arteriovenous malformation  
Hemorrhage  
Subdural hematoma  
Multiple sclerosis  
Atrophy of the brain

### Heart

Myocardial ischemia or infarction  
Ventricular dysfunction or enlargement  
Valvular disease  
Intracardiac thrombus  
Pericarditis or effusion  
Cardiac or pericardial masses  
Ventricular dilation or hypertrophy  
Congenital heart defects (e.g., septal defects)  
Diseases of the great vessels

### Other

Tumor (primary or metastatic)  
Abscess  
Edema



Bone destructive lesion  
Joint disorder  
Degenerative vertebral discs  
Kidney stones or gallstones

**notes**

**mammography** (Mammogram)

**Type of test** X-ray

**Normal findings**

Category 1: Negative

Category 2: Benign findings noted

Category 3: Probably benign findings; short-term follow-up is suggested

Category 4: Suspicious findings; further evaluation is indicated

Category 5: Cancer is highly suspected

Category 6: Known breast cancer

Category 0: Abnormality noted for which more imaging is recommended

**Test explanation and related physiology**


Mammography can identify cancers of the breast. In many cases, breast cancers can be detected before they become palpable. Radiographic signs of breast cancer include fine, stippled, clustered calcifications (white specks on the breast radiographs); a poorly defined, spiculated mass; asymmetric density; and skin thickening.

The detection rate for breast cancer with mammography is greater than 85%. This means that fewer than 15% of breast cancers are missed at mammography. Mammography also can detect other diseases of the breast such as acute suppurative mastitis, abscess, fibrocystic changes, cysts, benign tumors (e.g., fibroadenoma), and intraglandular lymph nodes.

Mammograms can be more thoroughly evaluated with magnified views, deeper views, or breast ultrasonography (see p. 181). Mammograms can be performed using analog radiographs or digital technology (*digital mammography*).

Mammography can also be used to locate a mammographically identified (i.e., not palpable) lesion for surgical or nonsurgical biopsy. Nonsurgical needle biopsy with a *stereotactic biopsy* device is the least invasive manner of obtaining tissue from a nonpalpable mammographic abnormality. For this procedure, the patient is placed prone on a specialized table ([Figure 29](#)). The mammogram is connected to a computer that can identify the exact location of the mammographic abnormality.

*Breast tomography (three-dimensional mammography)* through different thicknesses of the breast tissue increases sensitivity of the

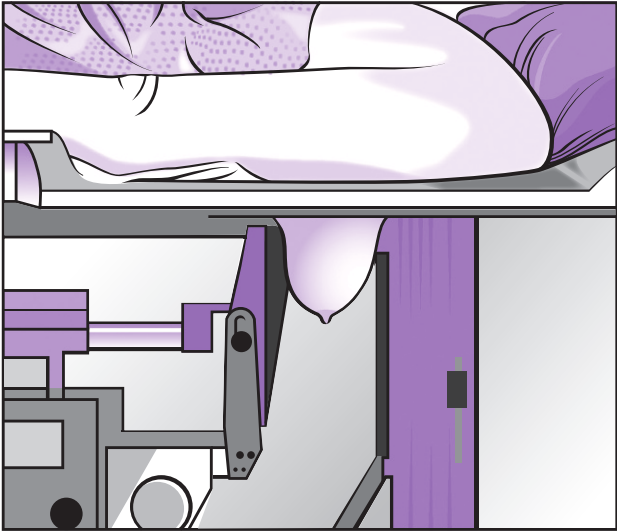

**FIG. 29** Stereotactic breast biopsy. The patient is positioned on the table with the breast pendulous through the aperture. The breast is compressed, with the target lesion centered in the biopsy window.

M

test. Unfortunately, this technique is too expensive for screening nonsymptomatic women.

The frequency and ages of asymptomatic women that benefit most from screening mammography is presently debated. See [Table 30](#).

### Contraindications

- Patients who are pregnant unless the benefits outweigh the risks of fetal damage
- Women younger than age 25 years

### Interfering factors

- Talc powder and deodorant can give the impression of calcification within the breast.
- Jewelry worn around the neck can preclude total visualization.
- Breast augmentation implants prevent total visualization of the breast. However, implants can be displaced so the native breast tissue can be imaged.
- Previous breast surgery can alter or distort the findings.

**TABLE 30** Summary of mammography screening guidelines

Age (yr)	USPSTF	ACS	ACOG	ACR	ASBS
< 50	2016 Consider mammograms.	2015 <i>Women aged 40-44yrs</i> should consider mammography. <i>Women aged 45-49yrs</i> should get annual mammograms.	2011 Yearly mammography	2010 Yearly mammography	2015 <i>Women aged</i> <i>40-44yrs</i> should consider mammography. <i>Women aged</i> <i>45-54yrs</i> should get yearly mammography.
50-75	Biennial screening mammography	<i>Women aged 50-54yrs</i> should get mammogram every year. <i>Women aged 55 yrs and</i> <i>older</i> should switch to mammography every 2 years or have the choice to continue yearly mammography.	Yearly mammography	Yearly mammography	Mammography yearly or every 2 yr for women aged 55 yrs and older

> 75	No screening	Mammography yearly or every 2 yrs as long as a woman is in good health and is expected to live 10 more years or longer	Women should, in consultation with their physicians, decide whether or not to continue mammographic screening.	Screening with mammography should stop when life expectancy is < 5-7 yrs on the basis of age or comorbid conditions.	Mammography every 2 yrs after the age of 75 yrs if the estimated life expectancy is > 10 yrs.
------	--------------	--	--	--	---

---

*ACOG*, American College of Obstetricians and Gynecologists; *ACR*, American College of Radiology; *ACS*, American Cancer Society; *ASBS*, American Society of Breast Surgeons; *USPSTF*, U. S. Preventative Services Task Force.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Inform the patient that some discomfort may be experienced during breast compression. This compression allows better visualization of the breast tissue. Assure the patient that the breast will not be harmed by the compression. Premenopausal women with very sensitive breasts may choose to schedule their mammogram 1 to 2 weeks after their menses to reduce any discomfort caused by compression.
- PT** Tell the patient that no fasting is required.
- PT** Explain to the patient that a minimal radiation dose will be used during the test.
- PT** Instruct the patient to disrobe above the waist and put on an x-ray gown.
  - Markers will be placed on any skin bump that may be interpreted as an abnormality on the x-ray image.

### During

- Note the following procedural steps:
  1. The procedure takes place in the radiology department or in a breast center with a mammogram machine.
  2. One breast is placed on the x-ray plate.
  3. The x-ray cone is brought down on top of the breast to compress it gently between the broadened cone and the x-ray plate.
  4. The x-ray image is obtained. This is the craniocaudal view.
  5. The x-ray plate is turned about 45 degrees medially and placed on the inner aspect of the breast.
  6. The broadened cone is brought in medially and again gently compresses the breast. This creates the mediolateral view.
  7. Occasionally other views, such as *direct lateral* (90 degree) or *magnified spot views*, are obtained to visualize more clearly an area of suspicion.
- Note that mammography is performed by a radiologic technologist in approximately 10 minutes. The x-ray images are interpreted by a radiologist.
- PT** Tell the patient that some discomfort may be caused by the pressure required to compress the breast tissue while the x-ray images are being taken.

**After**

**PT** Take this opportunity to instruct the patient in breast self-examination.

**Abnormal findings**

Breast cancer

Benign tumor (e.g., fibroadenoma)

Breast cyst

Fibrocystic changes

Breast abscess

Suppurative mastitis

**notes**

## maternal screen testing (Maternal triple screen, Maternal quadruple screen)

**Type of test** Blood or urine

### Normal findings

Low probability of fetal defects

### Test explanation and related physiology

These tests are provided to pregnant women early in pregnancy to identify potential birth defects or serious chromosomal or genetic abnormalities. These screening tests may indicate the potential for the presence of fetal defects (particularly trisomy 21 [Down syndrome] or trisomy 18). They may also indicate increased risk for neural tube defects (e.g., myelomeningocele, spina bifida) or abdominal wall defects (omphalocele or gastroschisis).

Maternal screening is routinely offered to all pregnant women, usually in their second trimester of pregnancy. Patients must understand that this is a screening test, not a diagnostic test. If the screening test results are positive, more accurate definitive testing such as blood tests, chorionic villus sampling (see p. 241), or amniocentesis (see p. 49) is recommended.

There are several variations of this test available:

- The *double test* measures two markers: human chorionic gonadotropin (hCG) (see p. 515) and alpha-fetoprotein (AFP; see p. 41).
- The *triple test* (maternal triple screen test) measures three markers: hCG, AFP, and estriol (see p. 386). AFP is produced in the yolk sac and fetal liver. Unconjugated estriol and hCG are produced by the placenta.
- The *quadruple test* measures four markers: hCG, AFP, estriol, and inhibin A.
- The *fully integrated screen test* measures AFP, estriol, fetal nuchal translucency (see p. 685), beta and total hCG, and pregnancy-associated plasma protein-A (PAPP-A; see p. 730).

The maternal triple screen test offers a 50% to 80% chance of detecting pregnancies with trisomy 21 compared with AFP alone, which has only a 30% chance of detection. The quadruple screen is now routinely recommended and is combined with fetal nuchal translucency (FTN) (see pelvic ultrasonography; see p. 685). These tests are most accurately performed during the second trimester of pregnancy, more specifically between the



14th and 24th weeks (ideal 16th and 18th weeks). The use of ultrasound to accurately indicate gestational age improves the sensitivity and specificity of maternal serum screening.

First trimester screening for genetic defects is an option for pregnant women. This testing would include fetal nuchal translucency (see pelvic ultrasonography; see p. 685), combined with the beta subunit of hCG (beta hCG; see p. 515), and pregnancy-associated plasma protein-A (PAPP-A; see p. 730). A low level of PAPP-A may indicate an increased risk for having a stillborn baby. These tests have detection rates comparable to standard second-trimester triple screening.

With trisomy 21, second trimester absolute maternal serum levels of AFP and unconjugated estriol are about 25% lower than normal levels, and maternal serum hCG is approximately two times higher than the normal hCG level. The results of the screening are expressed in *multiples of median (MoM)*. AFP and urinary estriol ( $E_3$ ) values during pregnancies with trisomy 21 are lower than those associated with normal pregnancies, which means that values lower than the mean are lower than 1 MoM. The hCG value for trisomy 21 is higher than 1 MoM. The MoM, fetal age, and maternal weight are used to calculate the possible risk for chromosomal abnormalities (e.g., trisomy 21). All of the previously named maternal screening tests are discussed elsewhere in this book. For the sake of thoroughness, inhibin A is discussed here.

*Inhibin A* levels in maternal serum remain relatively constant through the 15th to 18th week of pregnancy. Inhibin A is important in the control of fetal development. Maternal serum levels of inhibin A are twice as high in pregnancies affected by trisomy 21 as in unaffected pregnancies. The discovery of this fact led to the inclusion of inhibin A in the serum screening tests for trisomy 21. Inhibin A concentrations are significantly lower in women with normal pregnancies than in women with pregnancies that result in spontaneous abortions. Furthermore, circulating concentrations of inhibin A appear to reflect tumor mass for certain forms of ovarian cancer (see p. 194, tumor markers).

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient.

- Allow the patient to express her concerns and fears regarding the potential for birth defects.

## 614 maternal screen testing

### During

- Most of these tests can be done with a venous blood sample in a red-top tube. The hCG and estriol can also be tested by collecting a urine sample.

### After

- PT** Provide the results to the patient during a personal consultation.
- Allow the patient to express her concerns if the results are positive.
- PT** Assist the patient in scheduling and obtaining more accurate diagnostic testing if the results are positive.

### Abnormal findings

Positive screening tests (trisomy 21, trisomy 18, neural tube defects, abdominal wall defects)

### notes

## measles rubeola antibody

**Type of test** Blood

### Normal findings

Negative

### Test explanation and related physiology

The measles virus is an RNA paramyxovirus and is not the virus that causes the German measles; see rubella (p. 793). Although it is most usually a self-limiting disease, the virus can easily be spread (by respiratory droplets) to nonimmune pregnant women and cause preterm delivery or spontaneous abortion.

Testing for measles virus includes serologic identification of immunoglobulin G (IgG) and IgM antibodies. IgG elevation represents a previous infection or prior immunization. IgM indicates an acute infection or prior immunization. A fourfold rise in IgM indicates a current infection.

This test is used to diagnose measles in patients with a rash or viral syndrome when the diagnosis cannot be made clinically. Even more important, however, this test is used to establish and document immunity by previous measles infection or by previous vaccination. Populations commonly tested to document immunity include college students, health care workers, and pregnant women.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**PT** Inform the patient when to return for a follow-up rubeola titer if indicated.

**PT** If the results are negative for immunity, recommend immunization. For women of childbearing age, vaccination should precede future pregnancy.

### Abnormal findings

Active measles virus infection

Previous measles virus infection leading to immunity

notes

## Meckel diverticulum nuclear scan

**Type of test** Nuclear medicine

### Normal findings

No increased uptake of radionuclide in the right lower quadrant of the abdomen

### Test explanation and related physiology

Meckel diverticulum usually occurs in the ileum approximately 2 feet proximal to the ileocecal valve. Approximately 20% to 25% of Meckel diverticulum is lined internally by ectopic gastric mucosa. This gastric mucosa can secrete acid and cause ulceration of the intestinal mucosa nearby. Bleeding, inflammation, and intussusception are other potential complications of this congenital abnormality.

Both normal gastric mucosa within the stomach and ectopic gastric mucosa in Meckel diverticulum concentrate  $^{99m}\text{Tc}$  pertechnetate. One can then expect to see a hot spot in the right lower quadrant of the abdomen at about the same time as the normal stomach mucosa is visualized. This is a very sensitive and specific test for this congenital abnormality.

It is possible that Meckel diverticulum is present but contains no ectopic gastric mucosa within. Usually these are not symptomatic. No concentration of radionuclide will occur within the diverticulum. This test is not helpful in these cases.

Other conditions can simulate a hot spot compatible with Meckel diverticulum containing ectopic gastric mucosa. Usually these are associated with inflammatory processes within the abdomen (e.g., appendicitis or ectopic pregnancy).

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Advise the patient to refrain from eating or drinking anything for 6 to 12 hours before the examination.
- A histamine  $\text{H}_2$ -receptor antagonist is usually given for 1 to 2 days before the scan. This blocks secretion of the radionuclide from the ectopic gastric mucosa and improves visualization of Meckel diverticulum.

### During

- The patient lies in a supine position, and a large-view nuclear detector camera is placed over the patient's abdomen to identify concentration of nuclear material after intravenous injection.
- Images are taken at 5-minute intervals for 1 hour.
- Patients may be asked to lie on the left side to minimize the excretion of the radionuclide from the normal stomach, flooding the intestine with radionuclide and precluding visualization of Meckel diverticulum.
- Occasionally glucagon is provided to prolong intestinal transit time and avoid downstream contamination with the radionuclide.
- Occasionally gastrin is given to increase the uptake of the radionuclide by the ectopic gastric mucosa.
- There is no pain associated with this test.

### After

- The patient is asked to void, and a repeat image is obtained. This is to ensure that Meckel diverticulum has not been hidden by a distended bladder.

**PT** Because only tracer doses of radioisotopes are used, inform the patient that no precautions need to be taken by others against radiation.

### Abnormal findings

Meckel diverticula

notes

## mediastinoscopy

**Type of test** Endoscopy

### Normal findings

No abnormal mediastinal lymph node tissue

### Test explanation and related physiology

Mediastinoscopy is a surgical procedure in which a rigid mediastinoscope (a lighted instrument scope) is inserted through a small incision made at the suprasternal notch. The scope is passed into the superior mediastinum to inspect the mediastinal lymph nodes and to remove biopsy specimens. Because these lymph nodes receive lymphatic drainage from the lungs, their assessment can provide information on intrathoracic diseases (e.g., carcinoma, granulomatous infections, and sarcoidosis). Tumors occurring in the mediastinum (e.g., thymoma or lymphoma) can also be biopsied through the mediastinoscope.

### Potential complications

- Puncture of the esophagus, trachea, or blood vessels

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Ensure that the physician has obtained the informed consent for this procedure.
- Check whether the patient's blood needs to be typed and crossmatched.
- Provide preoperative care as with any other surgical procedure.
- Keep the patient NPO (nothing by mouth) after midnight on the day of the test.

#### During

- Note the following procedural steps:
  1. The patient is taken to the operating room for this surgical procedure.
  2. The patient is placed under general anesthesia.
  3. An incision is made in the suprasternal notch.
  4. The mediastinoscope is passed through this neck incision and into the superior mediastinum.
  5. The lymph nodes are biopsied.
  6. The scope is withdrawn, and the incision is sutured closed.

- Note that this procedure is performed by a surgeon in approximately 1 hour.
- PT** Inform the patient that he or she is asleep during the procedure.

### After

- Provide postoperative care as with any other surgical procedure.
- Assess the patient for mediastinal crepitus on auscultation, which may indicate mediastinal air from pneumothorax or the bronchus or esophagus.
- Note that distended neck veins and pulsus paradoxus may indicate lack of cardiac filling caused by a large mediastinal hematoma.
- Assess the patient for cough or shortness of breath, which may indicate a pneumothorax.
- Observe the patient for hypotension and tachycardia, which may indicate bleeding from the biopsy site or the great vessels.
- Evaluate the patient for hoarseness, which may indicate injury to the recurrent laryngeal nerve.
- Assess the patient for fever, chills, and sepsis, which may indicate mediastinitis from infection.

### Abnormal findings

Lung cancer  
 Metastasis  
 Sarcoidosis  
 Thymoma  
 Tuberculosis  
 Hodgkin disease  
 Lymphoma  
 Infection

### notes

**metal testing** (lead, arsenic, mercury)**Type of test** Blood**Normal findings**

&lt; 10 mcg/dL for lead

**Critical values for lead**Child ( $\leq 15$  years):  $\geq 20$  mcg/dLAdult ( $\geq 16$  years):  $\geq 70$  mcg/dL**Test explanation and related physiology**

Lead exposure, indicated by elevated blood lead levels, can result in permanent damage of almost all parts of the body. Children younger than 6 years of age are the most likely to be exposed and affected by lead. Blood lead levels are the best test for detecting and evaluating recent acute and chronic exposure. Blood lead samples are used to screen for exposure and to monitor the effectiveness of treatment.

Toxicity to heavy metals (e.g., *arsenic, mercury, cadmium, copper, cobalt*) and trace elements (e.g., *aluminum, chromium, antimony, bismuth, selenium, silver, zinc*) can also be tested. Suitable specimens for metal testing include blood, serum, urine, bones, teeth, tissue, or hair.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: verify with laboratory
- A finger stick can be performed to obtain 1 mL of blood.
- The blood sample is usually sent to a central diagnostic laboratory. The results are available in about 7 to 10 days.

**Abnormal findings****▲ Increased levels**

Exposure to lead, other heavy metals, and trace elements

**notes**



**metanephrine, plasma free** (Fractionated metanephrines)

**Type of test** Blood

**Normal findings**

Normetanephrine: < 0.5 nmol/L or 18-111 pg/mL by high-performance liquid chromatography (HPLC)

Metanephrine: < 0.9 nmol/L or 12-60 pg/mL by HPLC  
(Results vary among laboratories.)

**Test explanation and related physiology**

Pheochromocytomas, although rarely a cause of hypertension, are potentially lethal tumors. They produce several catecholamines that can cause episodic or persistent hypertension that is unresponsive to conventional treatment. The current diagnosis of pheochromocytoma depends on biochemical evidence of catecholamine overproduction by the tumor. This blood test measures the amount of metanephrine and normetanephrine, which are metabolites of epinephrine and norepinephrine, respectively.

The high sensitivity of plasma free metanephrine testing provides a high negative predictive value to the test. This means that if the concentrations of the free metanephrines are normal in the blood, then it is very unlikely that a patient has a pheochromocytoma. In about 80% of patients with pheochromocytoma, the magnitude of increase in plasma free metanephrines is so large that the tumor can be confirmed with close to 100% probability. Intermediate concentrations of normetanephrine and metanephrine are considered indeterminate. Urinary testing may clarify indeterminate findings.

When interpreting results, the following may be helpful:

- Any sample in which the concentrations of *both* normetanephrine and metanephrine are less than the upper reference range limit should be considered normal, and the presence of pheochromocytoma is highly unlikely.
- Any sample in which the concentrations of *either* normetanephrine or metanephrine exceed their respective upper reference range limits should be considered elevated.
- Whenever the normetanephrine or metanephrine concentration exceeds the indeterminate range, the presence of pheochromocytoma is highly probable and should be located via imaging techniques. Pheochromocytoma suppression and provocative testing (see p. 692) may assist in identifying this tumor.

## 622 metanephrine, plasma free

### Interfering factors

- Increased levels may be caused by caffeine or alcohol.
- Vigorous exercise, stress, and starvation may cause increased metanephrine levels.
- ▼ Drugs that may cause *increased* metanephrine levels include epinephrine-containing or norepinephrine-containing drugs, levodopa, lithium, and nitroglycerin.
- ▼ Acetaminophen can interfere with HPLC testing of metanephrines and should be avoided for 48 hours before testing.

### Procedure and patient care

#### Before

- ▶ Explain the procedure to the patient.
- ▶ Explain the dietary and medicinal restrictions.

#### During

- Identify and minimize factors contributing to patient stress.
- Physical and emotional stress may alter test results.
- The patient may be asked to lie down and rest quietly for 15 to 30 minutes before sample collection.
- The blood sample may be collected while supine.
- Collect a venous blood sample in a chilled lavender-top or pink-top tube. Invert to mix with preservatives.

#### After

- Apply pressure to the venipuncture site.
- Send the specimen to the laboratory as soon as the test is completed.

### Abnormal findings

#### ▲ Increased levels

Pheochromocytoma

notes

## methemoglobin (Hemoglobin M)

**Type of test** Blood

### Normal findings

0.06-0.24 g/dL or 9.3-37.2  $\mu\text{mol/L}$  (SI units)

0.4%-1.5% of total hemoglobin

### Possible critical values

> 40% of total hemoglobin

### Test explanation and related physiology

Methemoglobin is formed during the production of normal adult hemoglobin. If oxygenation of the iron component in the protohemoglobin occurs without subsequent reduction of the heme iron back to its  $\text{Fe}^{+2}$  form as exists in normal hemoglobin, excess methemoglobin accumulates. The oxidized iron form in methemoglobin is unable to combine with oxygen to carry the oxygen to the peripheral tissues. Therefore the oxyhemoglobin dissociation curve is shifted to the left, resulting in cyanosis and hypoxia.

Methemoglobinemia can be congenital or acquired. Hemoglobin M disease is a genetic defect that results in a group of abnormal hemoglobins that are methemoglobins. Another genetic mutation can cause a deficiency in reduced nicotinamide adenine dinucleotide (NADH) methemoglobin reductase enzyme, which is required to deoxygenate methemoglobin to normal adult hemoglobin. These forms of methemoglobinemia occur in infants, are usually severe, are not amenable to treatment, and are often fatal.

Acquired methemoglobinemia is a result of ingestion of nitrates (e.g., from well water) or such drugs as phenacetin, sulfonamides, isoniazid, local anesthetics, and some antibiotics. This form of the disease commonly occurs in older individuals and results in an acute crisis that is treated effectively with ascorbic acid or methylene blue. Methylene blue, however, is contraindicated in glucose-6-phosphate dehydrogenase (G6PD) deficiency.

### Interfering factors

- Tobacco use and carbon monoxide poisoning are associated with increased methemoglobin levels.
- ✦ Drugs that may cause *increased* levels include some antibiotics, isoniazid, local anesthetics, and sulfonamides.

## 624 methemoglobin

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: green
- Methemoglobin is very unstable. Place the specimen in an ice slush immediately after collection.
- Be prepared to provide oxygen support and close monitoring in the event the patient becomes increasingly hypoxic.

### Abnormal findings

Methemoglobinemia (hereditary or acquired)

notes

## methylated septin 9 DNA (mSEPT9, ColoVantage)

**Type of test** Blood

### Normal findings

0.0005-50 ng DNA

### Test explanation and related physiology

Because of the inconvenience and discomfort associated with routine screening (see colonoscopy, p. 258; stool for occult blood testing, p. 844) for colorectal cancer, many choose not to undergo testing and thereby preclude the opportunity for the early detection of an intestinal cancer. Recently a blood test for the detection of methylated septin 9 DNA (mSEPT9) has been developed that, when the result is positive, is very sensitive for the presence of a colorectal cancer.

A positive test result means there is an increased likelihood for the presence of a colorectal cancer or polyp. Individuals with positive test results are encouraged to undergo a diagnostic colonoscopy. Not all individuals with colorectal cancer will have a positive test result. Therefore individuals with a negative result should follow usual colorectal cancer screening guidelines.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

### Abnormal findings

- ▲ **Increased levels**
  - Colorectal cancer
  - Colorectal polyps

notes

**microalbumin (MA)****Type of test** Urine**Normal findings**

MA: &lt; 2 mg/dL

MA/creatinine ratio:

Males: &lt; 17 mg/g creatinine

Females: &lt; 25 mg/g creatinine

**Test explanation and related physiology**

Microalbuminuria refers to an albumin concentration in the urine that is greater than normal but not detectable with routine protein testing. Normally, only small amounts of albumin are filtered through the renal glomeruli, and that small quantity can be reabsorbed by the renal tubules. However, when the increased glomerular permeability of albumin overcomes tubular reabsorption capability, albumin is spilled in the urine. Preceding this stage of disease is a period of microalbuminuria that would normally go undetected. Therefore MA is an early indication of renal disease.

For a patient with diabetes, the amount of albumin in the urine is related to duration of the disease and the degree of glycemic control. MA is the earliest indicator for the development of diabetic complications (nephropathy, cardiovascular disease [CVD], and hypertension). MA can identify diabetic nephropathy 5 years before routine protein urine tests. Patients with diabetes with elevated MA have a 5- to 10-fold increase in the occurrence of CVD mortality, retinopathy, and end-stage kidney disease.

It is recommended that all patient with diabetes older than the age of 12 years be screened annually for MA. This can be done through a spot urine specimen by using a semiquantitative Micral Urine Test Strip. If MA is present, the test should be repeated two more times. If two of three MA urine test results are positive, a quantitative measurement using a 24-hour urine specimen should be performed.

The presence of MA in people without diabetes is an early indicator of lower life expectancy from CVD and hypertension. Nondiabetic nephropathies also may be associated with microalbuminuria. Life insurance underwriters are increasingly using MA testing to indicate life expectancy.

## Interfering factors

- Urinary tract infection, blood, or acid–base abnormalities can cause elevated MA levels.
- ✚ Oxytetracycline may interfere with test results.

## Procedure and patient care

- See inside front cover for Routine Urine Testing.
  - Ensure that the patient does not have any acute infection or urinary bleeding that could cause a false-positive result.
  - If the urine specimen contains vaginal discharge or bleeding, a clean-catch or midstream specimen will be needed.
  - Ensure that the urine sample is at room temperature for testing.
  - If using a Micral Urine Test Strip:
    1. Dip the test strip into the urine for 5 seconds.
    2. Allow the strip to dry for 1 minute.
    3. Compare the strip with the color scale on the label.
  - If a 24-hour urine collection is requested, the specimen should be refrigerated.
- PT** If the results are positive, inform the patient that the test should be repeated in 1 week.

## Abnormal findings

### ▲ Increased levels

Diabetes mellitus  
Hypertension  
Cardiovascular disease  
Nephropathy  
Urinary bleeding  
Hemoglobinuria  
Myoglobinuria

notes

**microglobulin** (Beta<sub>2</sub>-microglobulin [ $\beta_2$ m], Alpha-1-microglobulin, Retinol binding protein)

**Type of test** Blood, urine, fluid analysis

**Normal findings**

Beta<sub>2</sub>-microglobulin:

Blood: 0.70-1.80 mcg/mL

Urine:  $\leq$  300 mcg/L

Cerebrospinal fluid (CSF): 0-2.4 mg/L

Retinol binding protein:

Urine:  $<$  163 mcg/24 hours

Alpha-1-microglobulin:

Urine:  $<$  50 years:  $<$  13 mg/g creatinine

$\geq$  50 years:  $<$  20 mg/g creatinine

**Test explanation and related physiology**

$\beta_2$ m is increased in patients with malignancies (especially B-cell lymphoma, leukemia, or multiple myeloma; see p. 194, tumor markers), chronic infections, and chronic severe inflammatory diseases.

$\beta_2$ m, alpha-1-microglobulin, and retinol binding proteins pass freely through glomerular membranes and are almost completely reabsorbed by renal proximal tubule cells. Because of extensive tubular reabsorption, under normal conditions very little of these proteins appear in the final excreted urine. Therefore an increase in the urinary excretion of these proteins indicates proximal tubular disease or toxicity or impaired proximal tubular function. Therefore these proteins are helpful in differentiating among various types of renal disease. In patients with aminoglycoside toxicity, heavy metal nephrotoxicity, or tubular disease, protein urine levels are elevated. Excretion is increased 100 to 1000 times normal levels in cadmium-exposed workers. This test is used to monitor these workers.

$\beta_2$ m is particularly helpful in the differential diagnosis of renal disease. If blood and urine levels are obtained simultaneously, one can differentiate glomerular from tubular disease. In glomerular disease, because of poor glomerular filtration, blood levels are high, and urine levels are low. In tubular disease, because of poor tubular reabsorption, the blood levels are low, and urine levels are high. Blood levels increase early in kidney transplant rejection.



Increased CSF levels of  $\beta_2m$  indicate central nervous system (CNS) involvement with leukemia, lymphoma, HIV, or multiple sclerosis.

### Interfering factors

- Results could be affected by recent nuclear imaging when  $\beta_2m$  testing is performed by radioimmunoassay.
- $\beta_2m$  is unstable in acid urine.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient to minimize anxiety.

#### During

##### *Blood*

- Collect a venous blood sample in a red-top tube.

##### *Urine*

- PT** Instruct the patient to begin the 24-hour collection after voiding. Follow guidelines on inside front cover.
- A random urine level can also be used as a specimen when corrected for creatinine.

#### After

- Apply pressure to the venipuncture site.

### Abnormal findings

#### ▲ Increased urine levels

Renal tubule disease  
 Drug-induced renal toxicity  
 Heavy metal-induced renal disease  
 Lymphomas, leukemia, myeloma

#### ▲ Increased serum levels

Lymphomas, leukemia, myeloma  
 Glomerular renal disease  
 Renal transplant rejection  
 Viral infections, especially HIV and cytomegalovirus  
 Chronic inflammatory processes

### notes

**mononucleosis rapid test** (Mononuclear heterophil test, Heterophil antibody test, Monospot test)**Type of test** Blood**Normal findings**

Negative (&lt; 1:28 titer)

**Test explanation and related physiology**

The mononucleosis rapid test is performed to diagnose infectious mononucleosis (IM), a disease caused by the Epstein-Barr virus (EBV). Detectable levels of the IM heterophile antibody can usually be expected to occur between the sixth and tenth day after the onset of symptoms. The level usually increases through the second or third week of illness and, thereafter, can be expected to persist, gradually declining over a 12-month period. The IM heterophile antibody has been associated with several diseases other than IM. These include leukemia, Burkitt lymphoma, pancreatic carcinoma, viral hepatitis, cytomegalovirus infections, and others.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**Abnormal findings**

Infectious mononucleosis

Chronic EBV infection

Chronic fatigue syndrome

Burkitt lymphoma

Some forms of chronic hepatitis

**notes**

**Mycoplasma pneumoniae antibodies, IgG and IgM****Type of test** Blood**Normal findings**

IgG

- ≤ 0.9 (negative)
- 0.91-1.09 (equivocal)
- ≥ 1.1 (positive)

IgM

- ≤ 0.9 (negative)
- 0.91-1.09 (equivocal)
- ≥ 1.1 (positive)

IgM by IFA

Negative (reported as positive or negative)

**Test explanation**

Several diseases have been associated with the *M. pneumoniae* infection, including pharyngitis, tracheobronchitis, pneumonia, and inflammation of the tympanic membrane. *M. pneumoniae* accounts for approximately 20% of all cases of pneumonia. Classically, it causes a disease that has been described as primary atypical pneumonia. These infections may be associated with cold agglutinin syndrome (see p. 253).

Positive IgM results are consistent with acute infection, although there may be some cross-reactivity associated with other *Mycoplasma* infections. A single positive IgG result only indicates previous immunologic exposure.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Transport the specimen immediately to the laboratory. Avoid undue cooling of the specimen that may lead to agglutination.

**Abnormal findings***Mycoplasma* infection

notes

**myelography** (Myelogram)

**Type of test** X-ray with contrast dye

**Normal findings**

Normal spinal canal

**Test explanation and related physiology**

By placing radiopaque dye into the subarachnoid space of the spinal canal, the contents of the canal can be radiographically outlined. Cord tumors, meningeal tumors, metastatic spinal tumors, herniated intravertebral discs, and arthritic bone spurs can be readily detected by this study. These lesions appear as canal narrowing or as varying degrees of obstruction to the flow of the dye column within the canal. The entire canal (from lumbar to cervical areas) can be examined. This test is indicated in patients with severe back pain or localized neurologic signs that suggest the canal as the location of these injuries. Because this test is usually performed by lumbar puncture (LP; see p. 576), all the potential complications of that procedure exist.

A *water-soluble* contrast material is now used for myelography. This contrast does not need to be removed at the end of the procedure because it is water soluble and will be completely resorbed by the blood and excreted by the kidneys.

**Contraindications**

- Patients with multiple sclerosis because exacerbation may be precipitated by myelography
- Patients with increased intracranial pressure because LP may cause herniation of the brain
- Patients with infection near the LP site because this may precipitate bacterial meningitis
- Patients who are allergic to iodinated dye

**Potential complications**

- Headache
- Meningitis
- Herniation of the brain
- Seizures
- Allergic reaction to iodinated dye
- Hypoglycemia or acidosis may occur in patients who are taking metformin (Glucophage) and receive iodinated dye.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- Ensure that the physician has obtained written and informed consent for this procedure.
  - Assess the patient for allergies to iodinated contrast dye.
  - Ascertain whether the patient has recently taken phenothiazines, tricyclic antidepressants, CNS stimulants, or amphetamines if metrizamide will be used. These medications should be avoided, because they could decrease the seizure threshold.
  - Have the patient empty the bladder and bowel before myelography if possible.
- PT** Explain to the patient that he or she must lie very still during the procedure.
- Note that food and fluid restrictions vary according to the type of dye used. Check with the radiology department for specific restrictions.
- PT** Inform the patient that he or she will be tilted into an up-and-down position on the table so that the dye can properly fill the spinal canal and provide adequate visualization.

### During

- Note the following procedural steps:
  1. LP (see p. 576) is performed.
  2. Fifteen milliliters of CSF if withdrawn, and 15 mL or more of radiopaque dye is injected into the spinal canal.
  3. The patient is placed in the prone position on the tilt table with the head tilted down.
  4. Representative x-ray images are taken.
  5. After myelography is performed, the needle is removed and a dressing is applied.
- Note that this procedure is done by a radiologist in approximately 45 minutes.
- Keep in mind that patient response varies from mild discomfort to severe pain.

### After

- Note that nursing interventions after the procedure depend on the type of contrast agent used.
- PT** Usually place the patient on bed rest with the head slightly elevated for several hours afterward, as indicated. Position the patient as specifically ordered by the physician in consultation with the radiologist.

## 634 myelography

- Monitor the patient's vital signs and ability to void.
- PT Encourage the patient to drink fluids to enhance excretion of the dye and to hasten replacement of CSF.
- See p. xxi for appropriate interventions concerning the care of patients with iodine allergy.

### Abnormal findings

Cord tumor

Meningeal tumor

Metastatic spinal tumor

Meningioma

Cervical ankylosing spondylosis

Arthritic lumbar stenosis

Herniated intravertebral discs

Arthritic bone spurs

Neurofibroma

Avulsion of nerve roots

Cysts

Astrocytoma

### notes

## myoglobin

**Type of test** Blood

### Normal findings

< 90 mcg/L or < 90 mcg/L (SI units)

### Test explanation and related physiology

Myoglobin is an oxygen-binding protein found in cardiac and skeletal muscle. Measurement of myoglobin provides an early index of damage to the myocardium in myocardial infarction or reinfarction. Increased levels, which indicate cardiac muscle injury or death, occur in about 3 hours. Although this test is more sensitive than creatine phosphokinase (CPK) isoenzymes (see p. 297), it is not as specific. The benefit of myoglobin over CPK-MB is that it may become elevated earlier in some patients.

Disease or trauma of the skeletal muscle also causes elevations in myoglobin. Because myoglobin is excreted in the urine and is nephrotoxic, urine levels must be monitored in patients with high levels. Myoglobin may turn the urine red.

### Interfering factors

- Recent administration of radioactive substances
- Increased myoglobin levels can occur after IM injections.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Myocardial infarction  
Skeletal muscle inflammation  
(myositis)  
Malignant hyperthermia  
Muscular dystrophy  
Skeletal muscle ischemia  
Skeletal muscle trauma  
Rhabdomyolysis

#### ▼ Decreased levels

Polymyositis

notes

**natriuretic peptides** (Atrial natriuretic peptide [ANP], Brain natriuretic peptide [BNP], N-terminal fragment of pro-brain [B-type] natriuretic peptide [NT-pro-BNP], C-type natriuretic peptide [CNP], Ventricular natriuretic peptide, CHF peptides)

**Type of test** Blood

### Normal findings

ANP: 22-77 pg/mL or 22-77 ng/L (SI units)

BNP: < 100 pg/mL

CNP: Yet to be determined

NT-pro-BNP: < 300 pg/mL

### Critical values

BNP: > 400 pg/mL (Heart failure likely)

### Test explanation and related physiology

Natriuretic peptides (NPs) are used to identify and stratify patients with congestive heart failure (CHF). NPs are neuroendocrine peptides that oppose the activity of the renin-angiotensin system. There are three major NPs: ANP, BNP, and CNP. BNP has been implicated in the pathophysiology of hypertension, congestive heart failure (CHF), and atherosclerosis. BNP is released in response to atrial and ventricular stretch, respectively. BNP, in particular, correlates well to left ventricular pressures. As a result, BNP is a good marker for CHF. The higher the levels of BNP are, the more severe the CHF. This test is used in urgent care settings to aid in the differential diagnosis of shortness of breath (SOB). If BNP is elevated, the SOB is because of CHF. If BNP levels are normal, the SOB is pulmonary and not cardiac. This is particularly helpful in evaluating SOB in patients with cardiac and chronic lung disease.

Furthermore, BNP is a helpful prognosticator and is used in CHF risk stratification. CHF patients whose BNP levels do not rapidly return to normal with treatment experience a significantly higher risk for mortality in the ensuing months than do those whose BNP levels rapidly normalize with treatment. In early rejection of heart transplants, BNP levels can be elevated.

In some laboratories, BNP is measured as an *N-terminal fragment of pro-brain (B-type) natriuretic peptide (NT-pro-BNP)*. Screening diabetics for BNP elevation to determine the risk for cardiac diseases is helpful. BNP is also elevated in patients with prolonged systemic hypertension and in patients with acute myocardial infarction (MI).



## Interfering factors

- BNP levels are generally higher in healthy women than healthy men.
- BNP levels are higher in older patients.
- BNP levels are elevated for 1 month postoperatively in patients who have had cardiac surgery.
- There are several different methods of measuring BNP. Normal values vary whether or not the whole protein or a BNP fragment protein is measured.
- ✦ Natreacor (nesiritide), a recombination form of the endogenous human peptide used to treat CHF, will *increase* plasma BNP levels for several days.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

## Abnormal findings

### ▲ Increased levels

Congestive heart failure  
Myocardial infarction  
Systemic hypertension  
Cor pulmonale  
Heart transplant rejection

notes

**neuron-specific enolase (NSE)****Type of test** Blood**Normal findings**

&lt; 8.6 mcg/L

**Test explanation and related physiology**

NSE is present in neuroendocrine cells and amine precursor uptake and decarboxylation (APUD) cells. Elevated NSE concentrations are observed in patients with neuroblastoma, pancreatic islet cell carcinoma, medullary thyroid carcinoma, pheochromocytoma, and other neuroendocrine tumors.

NSE levels are frequently increased in patients with small cell lung cancer (SCLC) and infrequently in patients with non-SCLC (see p. 194, tumor markers). NSE can be used to monitor disease progression and management in SCLC. Levels of NSE can occasionally be elevated in benign disorders such as pneumonia and benign hepatobiliary diseases.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**Abnormal findings****▲ Increased levels**

Small cell lung cancer  
Neuroblastoma  
APUD tumors

**notes**

**neutrophil antibody screen** (Granulocyte antibodies, Polymorphonucleocyte antibodies [PMN ab], Antigranulocyte antibodies, Antineutrophil antibodies, Neutrophil antibodies, Leukoagglutinin)

**Type of test** Blood

### Normal findings

Negative for neutrophil antibodies

### Test explanation and related physiology

Neutrophil antibodies are antibodies directed toward white blood cells (WBCs). They develop during blood transfusions or transplacental bleeds and sometimes in patients with autoimmune disorders. Patients who experience a transfusion reaction despite complete compatibility testing before blood administration should have a neutrophil antibody screen to see if WBC incompatibility is the source of the reaction. This test is most commonly a part of posttransfusion antibody screening, which is a battery of testing performed if a transfusion reaction is suspected. It is also used in infants in the evaluation of unexplained neutropenia and in patients with suspected or known autoimmune disease.

Most commonly, in blood transfusion reactions, the recipient has antibodies to the donor WBCs and experiences a fever during transfusion. More severe, however, is the reaction when the donor plasma contains antibodies to the recipient's WBCs. This nonhemolytic reaction can lead to severe transfusion reactions, including acute pulmonary failure (*transfusion-related acute lung injury [TRALI]*) and multiorgan system failure.

### Interfering factors

- Recent administration of dextran
- Recent administration of intravenous (IV) contrast media
- Blood transfusion in the past 3 months

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or lavender

### Abnormal findings

TRALI

Alloimmune neonatal neutropenia

Autoimmune neutropenia

notes

**neutrophil gelatinase–associated lipocalin (NGAL, Lipocalin-2)****Type of test** Blood**Normal findings**

No rise in NGAL from baseline. (Results vary according to testing methods.)

**Test explanation and related physiology**

NGAL is a predictor for acute kidney injury (AKI), previously referred to as acute renal failure, and chronic kidney disease (CKD). There are no early markers for acute or chronic renal disease. Serum creatinine levels rise only after there has been significant renal impairment and injury. This is particularly important in patients who have serious nonrenal disease (e.g., heart surgery, renal transplant, sepsis). In these patients, severe AKI increases morbidity and mortality of hospitalized patients.

A marked elevation in NGAL indicates that renal injury has occurred. NGAL concentrations rise 48 hours before a rise in creatinine is noted. NGAL can be detected in both urine and blood within 2 hours of a renal insult.

NGAL can be measured in the urine, plasma, or serum samples with enzyme-linked immunosorbent assay (ELISA) test kits. Results are available in less than 1 hour in a standard laboratory with conventional ELISA equipment. This is particularly helpful in an intensive care environment. Normal values vary according to which laboratory method is used and the patient's baseline glomerular filtration rate (GFR). NGAL varies inversely with the GFR.

NGAL measurements are being used increasingly in a variety of clinical situations leading to AKI (e.g., during cardiac surgery, kidney transplantation, contrast nephropathy, and hemolytic uremic syndrome) and in the intensive care setting. They are also useful in conditions leading to CKD (e.g., lupus nephritis, glomerulonephritis, obstruction, dysplasia, polycystic kidney disease, IgA nephropathy, renal dysplasia, obstructive uropathy, and glomerular and cystic diseases).

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

- Collect urine specimens at the same time each day for consecutive days.
- Results are compared with previous day's testing.

### Abnormal findings

#### ▲ Increased levels

Primary or secondary renal disease

notes

## newborn metabolic screening

**Type of test** Blood

### Normal findings

Negative

### Possible critical values

Positive for any one of the tests

### Test explanation and related physiology

Newborn screening (NBS) is the practice of testing every newborn for certain harmful or potentially fatal disorders that are not otherwise apparent at birth. The Recommended Uniform Screening Panel (RUSP) is a list of disorders recommended for testing by the government. In most states, this testing is mandatory. NBS tests take place in the first 12 to 48 hours of life. If these diseases are not accurately diagnosed and treated, they can cause cognitive impairment, severe illness, and premature death in newborns. Many of these are metabolic disorders, often called *inborn errors of metabolism*. Disorders that may be detected through screening include endocrine and hematologic diseases.

Shortly after a child's birth, a sample of blood is obtained from a heel stick, and the blood is analyzed. The sample, called a *blood spot*, is tested at a reference laboratory. It is generally recommended that the sample be taken *after* the first 24 hours of life. Some tests such as the one for phenylketonuria (PKU) may not be as sensitive until the newborn has ingested an ample amount of the amino acid phenylalanine, which is a constituent of both human and cow milk, and after the postnatal thyroid surge has subsided. This is generally after about 2 days.

*Tandem mass spectrometry* can detect the blood components that are elevated in certain disorders, and it is capable of screening for more than 20 inherited metabolic disorders with a single test. The following disorders are typically included in newborn screening programs:

- **PKU:** An inherited disease, PKU is characterized by deficiency of the enzyme phenylalanine hydroxylase, which converts phenylalanine to tyrosine. Phenylalanine is an essential amino acid necessary for growth; however, any excess must be degraded by conversion to tyrosine. An infant with PKU lacks the ability to make this necessary conversion. Thus phenylalanine accumulates in the body and spills over into the urine. If the amount of phenylalanine is not restricted

in infants with PKU, progressive mental retardation results. A low-phenylalanine diet will need to be followed throughout childhood and adolescence and perhaps into adult life. (Incidence: 1 in 10,000 to 25,000.)

- *Congenital hypothyroidism*: Affected babies without treatment experience retarded growth and brain development. If the disorder is detected early, a baby can be treated with oral doses of thyroid hormone to permit normal development. (Incidence: 1 in 4000.)
- *Galactosemia*: Babies with galactosemia lack the enzyme that converts galactose into glucose, a sugar the body is able to use. As a result, milk and other dairy products must be eliminated from the diet. Otherwise, galactose can build up and cause blindness, severe mental retardation, growth deficiency, and even death. (Incidence: 1 in 60,000-80,000.) There are several less severe forms of galactosemia that may be detected by newborn screening. These may not require any intervention.
- *Sickle cell anemia*: Sickle cell disease is an inherited blood disease in which red blood cells stretch into abnormal *sickle* shapes (see p. 820). This can cause episodes of pain, damage to vital organs (e.g., the lungs and kidneys), and even death. Young children with sickle cell anemia are especially prone to certain dangerous bacterial infections. (Incidence: about 1 in 500 African American births and 1 in 1000 to 1400 Hispanic American births.)
- *Biotinidase deficiency*: Babies with this condition do not have enough biotinidase, an enzyme that recycles biotin (one of the B vitamins) in the body. This deficiency may cause seizures, poor muscle control, immune system impairment, hearing loss, mental retardation, coma, and even death. If the deficiency is detected early, however, problems can be prevented by biotin administration. (Incidence: 1 in 126,000.)
- *Congenital adrenal hyperplasia*: This is actually a group of disorders resulting in a deficiency of adrenal hormones. It can affect the development of the genitals and may cause death. Lifelong treatment through hormone supplementation manages the condition. (Incidence: 1 in 12,000.)
- *Maple syrup urine disease (MSUD)*: Babies with MSUD are missing an enzyme needed to process the amino acids leucine, isoleucine, and valine (present in protein-rich foods such as milk, meat, and eggs) that are essential for the body's normal growth. When these are not processed properly, they can build up in the body, causing urine to smell like maple syrup or sweet, burnt sugar. These babies usually have little

appetite and are extremely irritable. If not detected and treated early, MSUD can cause mental retardation, physical disability, and even death. A carefully controlled diet free of high-protein foods can prevent these outcomes. (Incidence: 1 in 250,000.)

- *Homocystinuria*: This metabolic disorder results from a deficiency in cystathionine  $\beta$ -synthase, responsible for the metabolism of methionine and homocysteine. If untreated, it can lead to dislocated lenses of the eyes, mental retardation, skeletal abnormalities, and hypercoagulability. However, a special diet combined with dietary supplements may help prevent most of these problems. (Incidence: 1 in 50,000-150,000.)
- *Tyrosinemia*: Babies with this disorder cannot metabolize tyrosine. If it accumulates in the body, it can cause mild retardation, language skill difficulties, liver problems, and even death from liver failure. A special diet and sometimes a liver transplant are needed to treat the condition. Early diagnosis and treatment seem to offset long-term problems. (Incidence: not yet determined.)
- *Cystic fibrosis*: This is an inherited disorder expressed in the lung and gastrointestinal tract that causes cells to release a thick mucus, leading to chronic respiratory disease, problems with digestion, and poor growth. There is no known cure; treatment involves trying to prevent the serious lung infections associated with it and providing adequate nutrition. (Incidence: 1 in 2000 white births.)
- *Toxoplasmosis*: Toxoplasmosis is a parasitic infection that can be transmitted through the mother's placenta to an unborn child. The disease-causing organism, which is found in undercooked meat, can invade the brain, eye, and muscle, possibly resulting in blindness and mental retardation. (Incidence: 1 in 1000.)

Certain other rare disorders can also be detected and may include Duchenne muscular dystrophy, HIV, and neuroblastoma. Hematologic disorders such as glucose-6-phosphate dehydrogenase (G6PD) deficiency and thalassemia can also be identified.

Most, but not all, states require newborns' hearing to be screened before they are discharged from the hospital. The hearing test involves placing a tiny earphone in the baby's ear and measuring his or her response to sound (see p. 393, evoked potentials).



## Interfering factors

- Premature infants may have false-positive results because of delayed development of liver enzymes.
- Infants tested before 24 hours of age may have false-negative results.
- Feeding problems may cause false-negative results.

## Procedure and patient care

### Before

**PT** Inform the parents about the purpose and method of the test.

- Assess the infant's feeding patterns before performing the test. An inadequate amount of protein ingested before performing the test can cause false-negative results.

### During

- Place a few drops of blood from a heel stick in each circle on the filter paper.

### After

**PT** Inform the parents that if test results are positive they will be notified by their healthcare provider, and further testing or treatment will be recommended.

## Abnormal findings

Metabolic diseases

Endocrine diseases

Hematologic diseases

## notes

**nicotine and metabolites** (Nicotine, Cotinine, 3-Hydroxy Cotinine, Nornicotine, Anabasine)

**Type of test** Urine/blood

**Normal findings**
**Urine**

	Unexposed nontobacco user (ng/mL)	Passive exposure (non-tobacco user) (ng/mL)	Abstinent user for > 2 weeks (ng/mL)	Active tobacco product user (ng/mL)
Nicotine	< 2	< 20	< 30	1000-5000
Cotinine	< 5	< 20	< 50	1000-8000
3-OH-Cotinine	< 50	< 50	< 120	3000-25,000
Nornicotine	< 2	< 2	< 2	30-900
Anabasine	< 3	< 3	< 3	3-500

**Serum**

	Unexposed nontobacco user (ng/mL)	Passive exposure (non-tobacco user) (ng/mL)	Abstinent user for > 2 weeks (ng/mL)	Active tobacco product user (ng/mL)
Nicotine	< 2	< 2	< 2	30-50
Cotinine	< 2	< 8	< 2	200-800
3-OH-Cotinine	< 2	< 2	< 2	100-500

**Test explanation and related physiology**

Nicotine is metabolized into cotinine and 3-hydroxy cotinine, which are measurable in urine and serum. The word “cotinine” is actually an anagram of “nicotine”—the eight letters are rearranged. In addition to nicotine and metabolites, tobacco products also contain other alkaloids (anabasine and nornicotine). These tests are used to assess compliance with smoking cessation programs and to qualify for surgical procedures. They are also used by insurance companies to determine whether the applicant

is a smoker. These tests can differentiate patient tobacco use as the following:

- Active user
- Abstinent > 2 weeks
- Passively exposed nonuser
- Unexposed nonuser

Cotinine and 3-hydroxy cotinine have *in vivo* half-lives of approximately 20 hours and are typically detectable for several days to up to 1 week after the use of tobacco. Because the level of these metabolites in the blood is proportionate to the amount of exposure to tobacco smoke, it is a valuable indicator of tobacco smoke exposure. Nicotine and its metabolites can be measured in the serum, urine, and other biofluids (most commonly the saliva). Cotinine is found in urine from 2 to 4 days after tobacco use. Serum or plasma testing is required when a valid urine specimen cannot be obtained or to detect recent use (within past 2 weeks). Blood cotinine will increase no matter how the tobacco is used (smoked, chew, dip, or snuff products). Nicotine levels have an *in vivo* half-life of approximately 2 hours, which is too short to be useful as a marker of smoking status.

Anabasine (only measured in the urine) is present in tobacco products but not nicotine replacement therapies. Nicotine, cotinine, 3-hydroxy cotinine, and nornicotine are also elevated by the use of any of the nicotine replacement gum, patch, or pill products. The presence of anabasine greater than 10 ng/mL or nornicotine greater than 30 ng/mL in urine indicates current tobacco use, irrespective of whether the subject is on nicotine replacement therapy. The presence of nornicotine without anabasine is consistent with use of nicotine replacement products. Heavy tobacco users who abstain from tobacco for 2 weeks exhibit urine nicotine values less than 30 ng/mL, cotinine values less than 50 ng/mL, anabasine values less than 2 ng/mL, and nornicotine values less than 2 ng/mL. Passive exposure to tobacco smoke can cause accumulation of nicotine metabolites in nontobacco users. Urine cotinine has been observed to accumulate up to 20 ng/mL from passive exposure. Neither anabasine nor nornicotine accumulates from passive exposure. Because hydration status and renal function may affect urinary cotinine results, a spot urine cotinine test is accompanied by a spot urine creatinine.

For smokers, another method of determining tobacco use is *expired carbon monoxide testing*. Again, a relatively short half-life (~4 hours) limits the reliability and accuracy. Furthermore, carbon monoxide testing is unable to detect the use of smokeless tobacco.

## Interfering factors

- Menthol cigarettes may increase cotinine levels because the menthol retains cotinine in the blood for a longer period of time.
- Diluted or adulterated urine may alter results.

## Procedure and patient care

### Before

- ▶ Explain the procedure to the patient and indicate the type of specimen needed.
- Obtain an accurate history of recent tobacco use.

### During

#### *Blood*

- Collect venous blood in a red-top, lavender-top (EDTA), or pink-top (K<sub>2</sub>EDTA) tube. See inside front cover for Routine Blood Testing.

#### *Urine*

- Obtain a random spot urine specimen of at least 5 mL. See inside front cover for Routine Urine Testing.
- Immediately transport the specimen to the laboratory.

#### *Saliva*

- Ask the patient to spit at least 1 mL of saliva into a spit container.
- Alternatively, dental gauze rolls can be placed in the mouth for 15 minutes and then placed in a storage container for transport.

### After

- Keep the specimens in a cool place if they cannot be transported to the laboratory immediately.

## Abnormal finding

Tobacco exposure

### notes

## 5'-nucleotidase

**Type of test** Blood

### Normal findings

0-1.6 units at 37° C or 0-1.6 units at 37° C (SI units)

### Test explanation and related physiology

5'-Nucleotidase is an enzyme specific to the liver. The 5'-nucleotidase level is elevated in patients with liver diseases, especially those associated with cholestasis. It provides information similar to alkaline phosphatase (ALP; see p. 29). However, ALP is not specific to the liver. When doubt as to the cause of an elevated ALP exists, 5'-nucleotidase is recommended. If that enzyme is elevated along with the ALP, the source of the pathology is certainly in the liver. If the 5'-nucleotidase is normal in the face of an elevated ALP, the source of pathology is outside the liver. Gamma-glutamyl transpeptidase (GGTP; see p. 435). is used similarly because it is also specific to the liver.

### Interfering factors

✘ Drugs that may cause *increased* 5'-nucleotidase levels include hepatotoxic agents.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Bile duct obstruction  
 Cholestasis  
 Hepatitis  
 Cirrhosis  
 Hepatic necrosis, ischemia, tumor  
 Hepatotoxic drugs

notes

## obstruction series

### Type of test

X-ray

### Normal findings

No evidence of bowel obstruction

No abnormal calcifications

No free air

### Test explanation and related physiology

The obstruction series is a group of x-ray images performed on the abdomen of patients with suspected bowel obstruction, paralytic ileus, perforated viscus, abdominal abscess, kidney stones, appendicitis, or foreign body ingestion. This series of images usually consists of at least two x-ray studies. The first is an *erect abdominal* image, which should include visualization of both diaphragms. The image is examined for evidence of free air under either diaphragm, which is pathognomonic for a perforated viscus. This view is also used to detect air–fluid levels within the intestine; the presence of an air–fluid level is compatible with bowel obstruction or paralytic ileus. Occasionally patients are too ill to stand erect. In this case an x-ray image can be taken with the patient in the left lateral decubitus position.

The second view in the obstruction series is usually a *supine abdominal* x-ray study. This x-ray is also called a *KUB* (*kidney, ureter, bladder*) film. A calcification within the course of the ureter could indicate a kidney or ureteral stone. An abdominal abscess may be seen as a cluster of tiny bubbles within one localized area. A small calcification in the right lower quadrant on the film of a patient complaining of pain in this quadrant may be an appendicolith. A gas-filled, distended bowel is compatible with bowel obstruction or paralytic ileus. The obstruction series can also be used to monitor the clinical course of patients with gastrointestinal (GI) disease.

Frequently, a *cross-table lateral* view of the abdomen is included in an obstruction series to detect abdominal aorta calcification, which often occurs in older patients.

Finally, the *supine abdominal* x-ray study can be used as a *scout image* before performing GI or abdominal x-ray studies that use contrast.

## Contraindications

- Patients who are pregnant unless the benefits outweigh the risk to the fetus

## Interfering factors

- Previous GI barium contrast study

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.

- Ensure that all radiopaque clothing has been removed.

**PT** Remind the patient that no GI contrast will be used.

### During

- Although the procedure varies from facility to facility, note that usually a supine abdominal x-ray image, erect abdominal image, and perhaps a lower erect chest image are taken. Often a cross-table lateral x-ray image is also included.
- Note that the obstruction series is performed in minutes in the radiology department by a radiologic technologist; however, it can be performed at the bedside with a portable x-ray machine. A radiologist interprets the images.

**PT** Tell the patient that no discomfort is associated with this study.

### After

- Note that no special aftercare is needed.

## Abnormal findings

Kidney stone

Bowel obstruction

Organomegaly

Presence of a foreign body

Bladder distention

Abdominal abscess

Perforated viscus

Abdominal aortic  
calcification

Appendicolithiasis

Paralytic ileus

Abdominal aortic aneurysm

Peritoneal effusion or ascites

Abnormal position of the  
kidneys

Soft tissue masses

## notes

**octreotide scan** (Carcinoid nuclear scan, Neuroendocrine nuclear scan)**Type of test**

Nuclear scan

**Normal findings**

No evidence of increased uptake throughout the body

**Test explanation and related physiology**

Octreotide scans are used to identify and localize neuroendocrine primary and metastatic tumors. These scans are indicated on patients with known neuroendocrine tumors.

Most neuroendocrine tumors have a somatostatin receptor on the cellular membrane. Octreotide is an analog of somatostatin. When combined with a radiopharmaceutical, the radiolabeled octreotide will attach to the somatostatin receptors of the neuroendocrine tumor cells. This test is used to identify primary and metastatic neuroendocrine tumors. It is also used to monitor the course of the disease.

Many different types of hormone-producing tumors can be detected by this scan, most notably carcinoid, gastrinoma, insulinoma, glucagonoma, pheochromocytoma, and small cell lung cancer. Other abnormalities can pick up octreotide, including granulomatous infections, rheumatoid arthritis, and non-hormonal cancers (breast, lymphoma, and non-small cell lung cancers).

**Contraindications**

- Patients who are pregnant or lactating because of risk of damage to the fetus or infant

**Interfering factors**

- Barium in the GI tract overlying the liver or spleen. Barium produces defects that may be mistaken for masses.

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Tell the patient that no fasting is required.
- PT** Assure the patient that he or she will not be exposed to large amounts of radiation because only tracer doses are used.
- If an iodinated radionuclide is to be used, ensure that the patient does not have an allergy to iodine.



- If an iodinated radionuclide is to be used, administer 5 drops of Lugol solution (iodine) daily for 3 days. This will avoid uptake of the radionuclide by the thyroid gland.
- PT** If the patient has been receiving octreotide as a form of anti-neoplastic treatment, this must be discontinued for 2 weeks before scanning. Inform the patient.

### During

- Note the following procedural steps:
  1. The patient is taken to nuclear medicine, where the radionuclide is administered intravenously.
  2. One hour after injection, a gamma ray detector or camera is successively placed over the entire body.
  3. The patient is placed in supine, lateral, and prone positions.
  4. The radionuclide image is recorded. Single-photon emission tomography (SPECT) images may also be performed.
  5. After 4 hours, the patient is given a strong laxative to clear the octreotide from the bowel.
  6. Repeat scanning is performed at 2, 4, 24, and 48 hours after administration of the octreotide.
- Note that the imaging procedure is performed by a trained technologist in approximately ½ hour. A physician trained in nuclear medicine interprets the results.
- PT** Inform the patient that the only discomfort associated with this procedure is the intravenous injection of the radionuclide.

### After

- PT** Because tracer doses of radioisotopes are used, inform the patient that no radiation precautions are needed.

### Abnormal findings

Carcinoid tumors

Neuroendocrine tumors

Granulomatous infections (e.g., sarcoidosis and tuberculosis)

### notes

## osmolality, blood (Serum osmolality)

### Type of test

Blood

### Normal findings

Adult/elderly: 285-295 mOsm/kg H<sub>2</sub>O or 285-295 mmol/kg (SI units)

Child: 275-290 mOsm/kg H<sub>2</sub>O

### Possible critical values

< 265 mOsm/kg H<sub>2</sub>O

> 320 mOsm/kg H<sub>2</sub>O

### Test explanation and related physiology

Osmolality measures the number of dissolved particles in serum/plasma per unit volume. As the amount of free water in the blood increases or the number of particles decreases per unit volume of serum, osmolality decreases. As the amount of water in the blood decreases or the number of particles per unit volume increases, osmolality increases. Osmolality increases with dehydration and decreases with overhydration.

The simultaneous use of urine osmolality (see p. 656) helps in the interpretation and evaluation of problems involving osmolality.

The serum osmolality test is useful in evaluating fluid and electrolyte imbalance. The test is very helpful in the evaluation of seizures, ascites, hydration status, acid-base balance, suspected antidiuretic hormone (ADH) abnormalities, and suspected poisoning. Osmolality is also helpful in identifying the presence of organic acids, sugars, and ethanol.

Osmolality can be predicted based on calculations of serum sodium, glucose, and blood urea nitrogen (BUN)—the three most important solutes in the blood. The equation is

$$\text{Osmolality} = 2 \times \text{Na} + \frac{\text{Glu}}{18} + \frac{\text{BUN}}{2.8}$$

The measured osmolality should not exceed the predicted by more than 10 mOsm/kg. A difference of more than 10 mOsm/kg is considered an *osmolal gap* or *delta gap*. Causes of a serum osmolal gap include mannitol, ethanol, methanol, ethylene glycol, and other toxins in very high concentration. Another measure providing similar data is the ratio of serum sodium to osmolality. Normally, the ratio of serum sodium, in mEq/L, to

serum osmolality, in mOsm/kg, is between 0.43 and 0.5. The ratio may be distorted in drug intoxication.

Osmolality may have a role in evaluation of coma patients. Values greater than 385 mOsm/kg  $H_2O$  are associated with stupor in patients with hyperglycemia. When values of 400 to 420 mOsm/kg are detected, grand mal seizures can occur. Values greater than 420 mOsm/kg can be lethal. The simultaneous use of urine osmolality helps in the interpretation and evaluation of problems with osmolality.

### Interfering factors

- Such diseases as cerebrovascular accident (stroke) or brain tumors may interfere with test results through inappropriate secretion of ADH.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- For pediatric patients, draw blood from a heel stick.

### Abnormal findings

#### ▲ Increased levels

Hyponatremia  
Dehydration  
Hyperglycemia  
Mannitol therapy  
Azotemia  
Uremia  
Ingestion of ethanol, methanol,  
or ethylene glycol  
Hyperosmolar nonketotic  
hyperglycemia  
Diabetes insipidus  
Hypercalcemia  
Renal tubular necrosis  
Severe pyelonephritis  
Ketosis  
Shock

#### ▼ Decreased levels

Hyponatremia  
Overhydration  
Syndrome of  
inappropriate  
antidiuretic  
hormone  
(SIADH) secretion  
Paraneoplastic  
syndromes  
associated with  
lung carcinoma

## osmolality, urine (Urine osmolality)

### Type of test

Urine

### Normal findings

12- to 14-hour fluid restriction: > 850 mOsm/kg H<sub>2</sub>O (SI units)  
 Random specimen: 50-1200 mOsm/kg H<sub>2</sub>O, depending on fluid intake, or 50-1200 mmol/kg (SI units)

### Test explanation and related physiology

Osmolality is the measurement of the number of dissolved particles in a solution. It is a more exact measurement of urine concentration than specific gravity.

Osmolality is used in the precise evaluation of the concentrating and diluting abilities of the kidney. With normal fluid intake and normal diet, a patient will produce a urine of about 500 to 850 mOsm/kg water. The normal kidney can concentrate a urine to 800 to 1400 mOsm/kg. With excess fluid intake, a minimal osmolality of 40 to 80 mOsm/kg can be obtained. With dehydration, the urine osmolality should be three to four times the plasma osmolality.

Osmolality is used in the evaluation of kidney function and the ability to excrete ammonium salts. Osmolality may be used as part of the urinalysis when the patient has glycosuria or proteinuria or has had tests that use radiopaque substances. In these situations, the *urine osmolar gap* increases because of other organic osmolar particles. The urine osmolar gap is the sum of all the particles predicted or calculated to be in the urine (electrolytes, urea, and glucose) compared with the actual measurement of the osmolality. The predicted or calculated urine osmolality can then be determined by urine levels of sodium, potassium, glucose, and urea nitrogen.

$$\text{Osmolality} = 2 \times ([\text{Na} + \text{K}]) + \frac{\text{Urea nitrogen}}{2.8} + \frac{\text{Glucose}}{18}$$

Normally, the osmolar gap is 80 to 100 mOsm/kg of H<sub>2</sub>O. The urine osmolality is more easily interpreted when the serum osmolality (see p. 654) is simultaneously performed. More information concerning the state of renal water handling or abnormalities of urine dilution or concentration can be obtained if urinary osmolality is compared with serum osmolality and if urine electrolyte studies are performed. Normally, the ratio of urine osmolality to serum osmolality is 1.0 to 3.0, reflecting a wide range of urine osmolality.

## Procedure and patient care

- See inside front cover for Routine Urine Testing.
- PT** Tell the patient that no special preparation is necessary for a random urine specimen.
- PT** Inform the patient that preparation for a fasting urine specimen may require a high-protein diet for 3 days before the test.
- PT** Instruct the patient to eat a dry supper the evening before the test and to drink no fluids until the test is completed the next morning.
- Preferably, collect a first-voided urine specimen for a random sample.
- Indicate on the laboratory slip the patient's fasting status.

## Abnormal findings

### ▲ Increased levels

SIADH secretion  
 Acidosis  
 Shock  
 Hyponatremia  
 Hepatic cirrhosis  
 Congestive heart failure  
 Addison disease

### ▼ Decreased levels

Diabetes insipidus  
 Hypercalcemia  
 Excess fluid intake  
 Renal tubular necrosis  
 Aldosteronism  
 Hypokalemia  
 Severe pyelonephritis

notes

**oximetry** (Pulse oximetry, Oxygen saturation)**Type of test**

Photodiagnostic

**Normal findings** $\geq 95\%$ **Possible critical values** $\leq 75\%$ **Test explanation and related physiology**

Oximetry is a noninvasive method of monitoring arterial blood oxygen saturation ( $\text{SaO}_2$ ). The  $\text{SaO}_2$  is the ratio of oxygenated hemoglobin to the total amount of hemoglobin. The  $\text{SaO}_2$  is expressed as a percentage; for example, a saturation of 95% indicates that 95% of the total hemoglobin attachments for oxygen have oxygen attached to them. The  $\text{SaO}_2$  is an accurate approximation of oxygen saturation obtained from an arterial blood gas study (see p. 106). By correlating the  $\text{SaO}_2$  and the patient's physiologic status, a close estimate of the partial oxygen pressure ( $\text{Po}_2$ ) can be obtained.


Oximetry is typically used for monitoring the patient's oxygenation status during the perioperative period (or any time of heavy sedation) and for patients receiving mechanical ventilation. This test is used in many clinical situations, such as pulmonary rehabilitation programs, stress testing, and sleep laboratories. Oximetry can be used to assess the body's response to various drugs, such as theophylline (which causes bronchodilation) and methacholine (which evokes bronchospasm in people with asthma). This test is commonly used to titrate levels of oxygen on hospitalized patients.

$\text{O}_2$  levels can also be measured in various body tissues. For example, monitors that continuously measure tissue  $\text{O}_2$  partial pressures can be attached to a small catheter placed in the brain, heart, or peripheral muscle. *Brain tissue oxygen testing and monitoring* is the most common use of this technology. It is used to monitor the condition of the brain after severe head trauma by measuring cerebral blood flow and pulmonary oxygenation. It is more accurate than intracranial pressure for indicating brain injury.

**Procedure and patient care****Before**

**PT** Explain the procedure to the patient.

**PT** Tell the patient that no fasting is required.

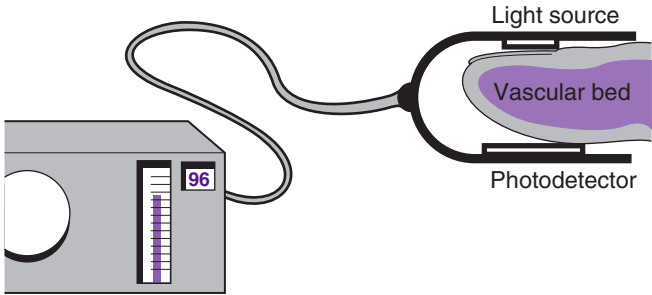

**FIG. 30** Oximetry. The pulse oximeter passes a beam of light through the tissue. The amount of light absorbed by the oxygen-saturated hemoglobin is measured by the sensor.

### During

- Rub the patient's earlobe, pinna (upper part of the ear), or fingertip to increase blood flow.
- Clip the monitoring probe or sensor to the ear or finger. A beam of light passes through the tissue, and the sensor measures the amount of light the tissue absorbs (Figure 30).
- Note that this study is usually performed by a respiratory therapist or nurse at the patient's bedside in a few seconds.

**PT** Tell the patient that no discomfort is associated with this study.

### After

- Note that no special aftercare is needed.

### Abnormal findings

#### ▲ Increased levels

Increased inspired  $O_2$   
Hyperventilation

#### ▼ Decreased levels

Inadequate  $O_2$  in inspired air  
Hypoxic lung diseases  
Hypoxic cardiac diseases  
Severe hypoventilation states

notes

**pancreatic elastase (PE, EL1, Elastase-1, Fecal Elastase)****Type of test** Stool**Normal findings**

Stool:

- > 200 mcg Elastase/g stool = Normal
- 100-200 mcg Elastase/g stool = Moderate to slight exocrine pancreatic insufficiency
- < 100 mcg Elastase/g stool = Severe exocrine pancreatic insufficiency

Serum:

- Less than 3.5 ng/mL

**Test explanation and related physiology**

PE is a digestive enzyme secreted by the pancreas. In this test, PE is used as a fecal marker to assess exocrine pancreatic function. This test is used to identify maldigestion caused by pancreatic insufficiency. The insufficiency can be caused by chronic pancreatitis, cystic fibrosis, or pancreatic resection. If possible, pancreatic enzyme supplementation should be discontinued before testing. EL1 can also be measured in the serum. Elevated serum levels may indicate pancreatitis. Hypothyroidism is associated with decreased levels.

**Interfering factors**

- ✱ Pancreatic enzyme supplementation can increase results.

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient. No preparation is needed.

**During**

- PT** Instruct the patient to transfer 5 g of stool to a sterile, plastic, screw-cap container provided by the laboratory. No preservative is needed.
- PT** Inform the patient that watery, diarrheal stool is not recommended for this test because some dilutions during the assay may give lower EL1 concentrations than are actually present.

**After**

- PT** Tell the patient results may take several days if the specimen is sent to an outside reference lab.



**Abnormal findings****▼ Decreased stool levels**

Pancreatic insufficiency (e.g.,  
chronic pancreatitis)

Pancreatic resection

Cystic fibrosis

Diabetes mellitus

Cholelithiasis

Osteoporosis

**▲ Increased serum levels**

Pancreatitis

**▼ Decreased serum levels**

Hypothyroidism

**notes**

## pancreatic enzymes (Pancreatic secretory test, Amylase, Lipase, Trypsin, Chymotrypsin)

**Type of test** Fluid analysis

### Normal findings

Volume: 2-4 mL/kg body weight

HCO<sub>3</sub> (bicarbonate): 90-130 mEq/L

Amylase: 6.6-35.2 units/kg

Trypsin-like immunoreactivity: 10-57 ng/mL

Trypsin:  $\geq$  1:96

Chymotrypsin: by report

### Test explanation and related physiology

Cystic fibrosis (CF) is an inherited disease characterized by abnormal secretion by exocrine glands within the bronchi, small intestines, pancreatic ducts, bile ducts, and skin (sweat glands). Because of this abnormal exocrine secretion, children with CF develop mucus plugs that obstruct their pancreatic ducts and can lead to significant malabsorption, steatorrhea, and diarrhea. The pancreatic enzymes (e.g., amylase, lipase, trypsin, and chymotrypsin) cannot be expelled into the duodenum and therefore are either completely absent or present only in diminished quantities within the duodenal aspirate. For the same reasons, bicarbonate and other neutralizing fluids cannot be secreted from the pancreas.

In this test, secretin and pancreozymin are used to stimulate pancreatic secretion of these enzymes and bicarbonate into the duodenum. The duodenal contents are then aspirated and examined for pH, bicarbonate, and pancreatic enzyme levels. Amylase is the most frequently measured enzyme.

*Trypsinogen*, another pancreatic exocrine enzyme, is measured in the serum as *trypsin-like immunoreactivity*. This test is used to support the diagnosis of chronic pancreatitis. Levels diminish as pancreatic exocrine function becomes increasingly impaired.

When any of these pancreatic enzymes are measured in the serum, they can reflect acute inflammation of the pancreas. Likewise, in patients with burned-out chronic pancreatitis, serum measurements of these pancreatic enzymes are low. Pancreatic enzymes can also be measured in the stool and are more easily accessed for testing. Most commonly *stool for trypsin and chymotrypsin* are tested and are a measure of pancreatic function.

Stool for elastase (see p. 660) also is a measure of pancreatic function.

A physician obtains the duodenal contents in approximately 2 hours in the x-ray department. Discomfort and gagging may occur during placement of the Dreiling tube.

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient and parents.

**PT** Instruct the adult patient to fast for 12 hours before testing.

- Determine pediatric fasting times according to the patient's age.

### During

- Note the following procedural steps:
  1. With the use of fluoroscopy, a Dreiling tube is passed through the patient's nose and into the stomach.
  2. The distal lumen of the tube is placed within the duodenum.
  3. The proximal lumen of the tube is placed within the stomach.
  4. Both lumens are aspirated. The gastric lumen is continually aspirated to avoid contamination of the gastric contents in the duodenum aspirate.
  5. A control specimen of the duodenal juices is collected for 20 minutes.
  6. The patient is tested for sensitivity to secretin and pancreozymin by low-dose intradermal injection.
  7. If no sensitivity is present, these hormones are administered intravenously. Secretin can be expected to stimulate pancreatic fluid and bicarbonate secretion. Pancreozymin can be expected to stimulate pancreatic enzyme secretion.
  8. Four duodenal aspirates are collected at 20-minute intervals and placed in the specimen container.
  9. Each specimen is analyzed for pH, volume, bicarbonate, and amylase levels.
- Note that a physician performs this test in the laboratory or at the patient's bedside.

### After

- Place the aspirated specimens on ice. Send them to the chemistry laboratory as soon as the test is completed.
- Remove the Dreiling tube after completion of the test. Give appropriate nose and mouth care.
- Allow the patient to resume a normal diet.

### Abnormal findings

▲ Increased levels

Acute pancreatitis

▼ Decreased levels

Cystic fibrosis

Sprue

Chronic pancreatitis

notes

## pancreatobiliary fish testing

**Type of test** Microscopic examination

### Normal findings

No chromosomal ploidy abnormalities

### Test explanation and related physiology

It is sometimes difficult to differentiate benign bile duct strictures from early pancreatobiliary cancer. When a stricture is identified on an endoscopic retrograde cholangiopancreatography (ERCP; see p. 364), cancer must be considered as a possible cause. If an obvious cancer is not seen at the time of ERCP, a brush is repeatedly swept along the bile duct to obtain duct surface cells for conventional cytology to identify cancer cells. In conventional cytology, the brushing specimens are placed on a slide and stained with a Papanicolaou (PAP) stain. Slides are then interpreted by a cytopathologist to determine whether they show features that are positive for malignancy, suspicious for malignancy, atypical (meaning there are cells that are not normal but cannot be definitely ascribed to a neoplastic process), or negative for malignancy.

With the use of fluorescence *in situ* hybridization (FISH) testing, three chromosome enumeration probes and a gene-specific probe to P16 tumor suppressor gene are able to determine whether more than one pair of chromosomes or P16 genes exists in the cells obtained from the brushings of the bile duct during ERCP. If extra copies of two or more of the chromosomes or P16 genes are evident, the cells are considered to be *polysomic*, which indicates a high chance of malignancy.

### Contraindications

- See ERCP.

### Potential complications

- See ERCP.

### Interfering factors

- Errors in obtaining a good specimen can influence results.
- Cytologic examination is always affected by physician interpretation.

## 666 pancreatobiliary fish testing

### Procedure and patient care

#### Before

- ▶ Explain the procedure to the patient.
- Obtain informed consent from the patient.
- Keep the patient NPO (nothing by mouth) as of midnight the day of the test.
- Follow the procedure for ERCP.

#### During

- During ERCP, a rounded brush is placed through the accessory lumen of the endoscope and passed repeatedly through the stricture.
- The brush is then swished in a cytology solution for FISH or directly smeared on a slide and preserved for conventional cytology.

#### After

- Follow the guidelines for ERCP.

### Abnormal findings

Sclerosing cholangitis

Biliary sclerosis

Tumor or strictures of the pancreatobiliary duct

notes

## Papanicolaou smear (Pap smear, Pap test, Cytologic test for cancer, Liquid-based cervical cytology [LBCC], ThinPrep)

**Type of test** Microscopic examination

### Normal findings

No abnormal or atypical cells

### Test explanation and related physiology

A Pap test can detect neoplastic cells in cervical and vaginal secretions. This test is based on the fact that normal cells and abnormal cervical and endometrial neoplastic cells are shed into the cervical and vaginal secretions. By examining these secretions microscopically, one can detect early cellular changes associated with premalignant conditions or an existing malignant condition. The Pap test is 95% accurate in detecting cervical carcinoma; however, its accuracy in the detection of endometrial carcinoma is only approximately 40%.

The Bethesda System includes evaluation of the following five components:

- Adequacy of specimen
- General categorization (optional)
- Interpretation or result
- Automated review and ancillary testing (e.g., human papillomavirus [HPV]; see p. 520)
- *Educational notes (optional)*
  - Negative for intraepithelial lesion or malignancy
  - Epithelial cell abnormalities

A common method of Pap test specimen collection is *liquid-based cervical cytology (LBCC)*. With this technique, the specimen obtained from the cervix is placed into a preservative solution instead of smearing it onto a slide as is done during conventional Pap smear testing (CPT). Any blood cells and debris are then isolated by centrifuge, leaving only cervical cells. A thin film of the residuum is then placed on a slide to be evaluated. The specimen can be split into two parts. The first is evaluated for cytopathology. In the event that cytologic abnormalities of undetermined significance are found that could be better elucidated with further testing, the cells in the second split specimen are used for that testing (to avoid having to obtain another cervical sample). For example, if cellular changes are found that may be related to HPV, the second split specimen is tested by real-time polymerase chain reaction (PCR) for HPV DNA (see p. 520).

## 668 papanicolaou smear

A slightly different technique called the *PapSpin* uses a special brush placed in a collection device and centrifuged to provide a cellular concentrate for microscopic examination.

U.S. Preventive Services Task Force (USPSTF) recommended that

- Women aged 21 to 65 years should get Pap tests no more than every 3 years.
- Women aged 30 to 65 years may extend the interval between screens to 5 years if they use HPV tests in conjunction with the Pap test.
- Women younger than 21 years should not be screened for cervical cancer, regardless of sexual history.
- Women older than 65 years should not be screened as long as they have had consistently normal Pap test results and are not at high risk for cervical cancer.

The guidelines apply to healthy women who do not have abnormal Pap test results. They do not apply to women who have a history of cervical cancer.

### Contraindications

- Patients currently having routine, normal menses because this can alter test interpretation
- Patients with vaginal infections: Cellular changes that may be misinterpreted as dysplastic may transiently occur during these infections.

### Interfering factors

- A delay in fixing a specimen allows the cells to dry, destroys the effectiveness of the stain, and makes interpretation difficult.
- Using lubricating jelly on the speculum can alter the specimen.
- Douching and tub bathing may wash away cellular deposits.
- Menstrual flow may alter test results.
- Infections may interfere with hormonal cytology.
- ✚ Such drugs as digitalis and tetracycline may alter the test results by affecting the squamous epithelium.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient not to douche or tub bathe during the 24 hours before the Pap test. (Some physicians prefer that patients refrain from sexual intercourse for 24-48 hours.)
- PT** Instruct the patient to empty her bladder.
- PT** Tell the patient that no fasting or sedation is required.




## During

- Note the following procedural steps:
  1. The patient is placed in the lithotomy position.
  2. A vaginal speculum is inserted to expose the cervix.
  3. Material is collected from the cervical canal by rotating a moist saline cotton swab or spatula within the cervical canal and in the squamocolumnar junction (Figure 31).
  4. The cells are immediately wiped across a clean glass slide and fixed either by immersing the slide in equal parts of 95% alcohol and ether or by using a commercial spray (e.g., Aqua Net hair spray). The secretions must be fixed before drying because drying will distort the cells and make interpretation difficult.
  5. If *liquid-based cervical cytology* is performed, the cervical specimen is placed in the fixative preservative solution. After being placed in this solution, cells can be evaluated any time within the next 3 weeks if kept frozen.
  6. The slide is labeled with the patient's name, age, and parity and with the date of her last menstrual period.
- Note that a Pap smear is obtained by a nurse or a physician in approximately 10 minutes.
- PT** Tell the patient that no discomfort, except for insertion of the speculum, is associated with this procedure.

## After

- PT** Inform the patient that usually she will not be notified unless further evaluation is necessary.

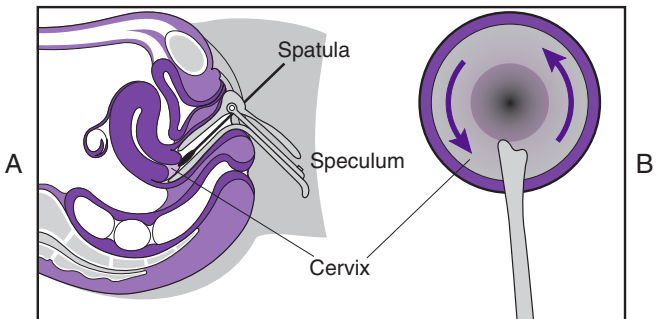

**FIG. 31** Papanicolaou (Pap) smear. **A**, The vaginal speculum is shown in position to allow direct visualization of the cervix. **B**, The cervix is scraped with the bifid end of a wooden spatula.

**670** papanicolaou smear

**Abnormal findings**

Cancer

Infertility

Sexually transmitted disease

HPV infection

Reactive inflammatory changes

Fungal infection

Parasitic infection

Herpes infection

**notes**

## paracentesis (Peritoneal fluid analysis, Ascitic fluid cytology, Peritoneal tap)

**Type of test** Fluid analysis

### Normal findings

Gross appearance: clear, serous, light yellow, < 50 mL

Red blood cells (RBCs): none

White blood cells (WBCs): < 300/ $\mu$ L

Protein: < 4.1 g/dL

Glucose: 70-100 mg/dL

Amylase: 138-404 units/L

Ammonia: < 50 mcg/dL

Alkaline phosphatase

Adult male: 90-240 units/L

Female < 45 years: 76-196 units/L

Female > 45 years: 87-250 units/L

Lactic dehydrogenase (LDH): similar to serum LDH

Cytology: no malignant cells

Bacteria: none

Fungi: none

Carcinoembryonic antigen (CEA): negative


### Test explanation and related physiology

Paracentesis is an invasive procedure entailing the insertion of a needle or catheter into the peritoneal cavity (Figure 32) for removal of ascitic fluid for diagnostic and therapeutic purposes.

*Diagnostically*, paracentesis is performed to obtain and analyze fluid to determine the etiology of the peritoneal effusion. Peritoneal fluid is classified as to whether it is a transudate or exudate. This is an important differentiation and is very helpful in determining the etiology of the effusion. *Transudates* are most frequently caused by congestive heart failure, cirrhosis, nephrotic syndrome, myxedema, peritoneal dialysis, and hypoproteinemia. *Exudates* are most often found in infectious or neoplastic conditions. However, collagen vascular disease, gastrointestinal (GI) diseases, trauma, and drug hypersensitivity also may cause an exudative effusion.

*Therapeutically*, this procedure is done to remove large amounts of fluid from the abdominal cavity.

The peritoneal fluid is usually evaluated for gross appearance, RBCs, WBCs, protein, glucose, amylase, ammonia, alkaline phosphatase, LDH, cytology, bacteria, fungi, and other tests

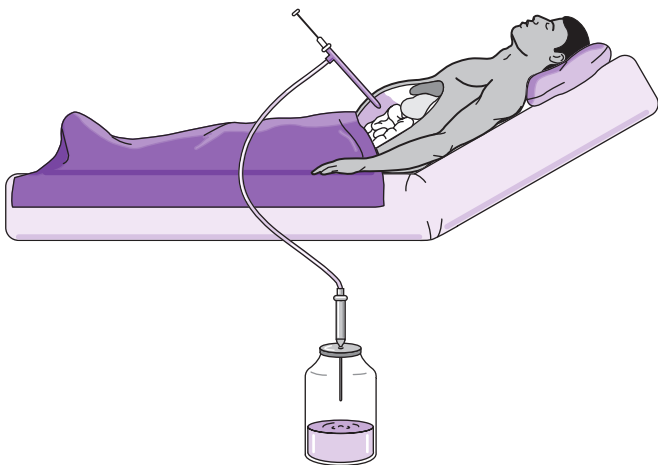

**FIG. 32** Paracentesis. A catheter is placed through the skin and abdominal muscle wall and into the peritoneal cavity containing free fluid.

(e.g., CEA levels). Urea and creatinine may be measured if there is a question that the fluid may represent urine from a perforated bladder.

### Contraindications

- Patients with coagulation or bleeding abnormalities
- Patients with only a small amount of fluid and extensive previous abdominal surgery

### Potential complications

- Hypovolemia if a large volume of peritoneal fluid is removed
- Peritonitis

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- Obtain informed consent for this procedure.
- PT** Tell the patient that no fasting or sedation is necessary.
- PT** Have the patient empty the bladder before the test.
- Measure abdominal girth.
- Obtain the patient's weight.
- Obtain baseline vital signs.

## During

- Note the following procedural steps:
    1. Place the patient in a high Fowler position in bed.
    2. Paracentesis is performed under a strict sterile technique.
    3. The needle insertion site is aseptically cleansed and anesthetized locally.
    4. A scalpel may be used to make an incision in the skin to allow the cannula or needle to enter.
    5. A trocar, cannula, or needle is threaded through the incision.
    6. Tubing is attached to the cannula. The other end of the tubing is placed in the collection receptacle (usually a container with a pressurized vacuum).
  - Note that this procedure is performed by a physician at the patient's bedside, in a procedure room, or in the physician's office in less than 30 minutes. Usually the volume removed is limited to about 4 L at any one time to avoid hypovolemia if the ascites is rapidly reaccumulated.
- PT** Although local anesthetics eliminate pain at the insertion site, tell the patient that he or she will feel a pressure-like discomfort as the needle is inserted.

## After

- All tests on peritoneal fluid should be performed immediately to avoid false results because of chemical or cellular deterioration.
- Place a small bandage over the needle site.
- Observe the puncture site for bleeding, continued drainage, or signs of inflammation.
- Measure the abdominal girth and weight of the patient; compare with baseline values.
- Monitor vital signs for evidence of hemodynamic changes.
- Because of the high protein content of ascitic fluid, albumin infusions may be ordered after paracentesis to compensate for protein loss.
- Occasionally ascitic fluid will continue to leak out of the needle tract after removal of the needle. A suture can stop that. If unsuccessful, a collection bag should be used.

## Abnormal findings

### Exudate

Lymphoma  
Carcinoma  
Tuberculosis

### Transudate

Hepatic cirrhosis  
Portal hypertension  
Nephrotic syndrome

## 674 paracentesis

### Exudate

Peritonitis

Pancreatitis

Ruptured viscus

### Transudate

Hypoproteinemia

Congestive heart failure

Abdominal trauma

Peritoneal bleeding

notes

## parathyroid hormone (PTH, Parathormone)

**Type of test** Blood

### Normal findings

Intact (whole): 10-65 pg/mL or 10-65 ng/L (SI units)

N terminal: 8-24 pg/mL

C terminal: 50-330 pg/mL

### Test explanation and related physiology

This test is useful in establishing a diagnosis of hyperparathyroidism and distinguishing nonparathyroid from parathyroid causes of hypercalcemia. Increased PTH levels are seen in patients with hyperparathyroidism (primary, secondary, or tertiary); in patients with nonparathyroid, ectopic PTH-producing tumors (pseudohyperparathyroidism); or as a normal compensatory response to hypocalcemia in patients with malabsorption or vitamin D deficiency.

It is important to measure serum calcium simultaneous with the measurement of PTH. Most laboratories have a PTH/calcium nomogram already made up, indicating what PTH level is considered normal for each calcium level.

Decreased PTH levels are seen in patients with hypoparathyroidism or as a compensatory response to hypercalcemia in patients with metastatic bone tumors, sarcoidosis, vitamin D intoxication, or milk-alkali syndrome. Of course, surgical ablation of the parathyroids is another cause of hypoparathyroidism.

### Interfering factors

- Recent injection of radioisotopes
- ✚ Drugs that *increase* PTH include anticonvulsants, isoniazid, lithium, rifampin, and steroids
- ✚ Drugs that *decrease* PTH include cimetidine and propranolol

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red
- Obtain a morning blood specimen, because diurnal rhythm affects PTH levels. (Check with the laboratory if the patient works at night.)
- Note that some laboratories require blood in an iced plastic syringe.

## 676 parathyroid hormone

- Obtain a serum calcium level determination at the same time if ordered. Serum PTH and serum calcium levels are important in the differential diagnosis.
- Indicate the time the blood was drawn on the laboratory slip because a diurnal rhythm affects test results.

### Abnormal findings

#### ▲ Increased levels

Hyperparathyroidism  
secondary to adenoma  
or carcinoma of the  
parathyroid gland  
Non-PTH-producing  
tumors (paraneoplastic  
syndrome)  
Lung carcinoma  
Kidney carcinoma  
Hypocalcemia  
Chronic renal failure  
Malabsorption syndrome  
Vitamin D deficiency  
Rickets  
Osteomalacia  
Congenital renal defect

#### ▼ Decreased levels

Hypoparathyroidism  
Hypercalcemia  
Metastatic bone tumor  
Sarcoidosis  
Autoimmune destruction of  
the parathyroid glands  
Vitamin D intoxication  
Milk-alkali syndrome  
Graves disease  
Hypomagnesemia  
DiGeorge syndrome

### notes



**parathyroid scan** (Parathyroid scintigraphy)

**Type of test** Nuclear medicine

**Normal findings**

No increased parathyroid uptake

**Test explanation and related physiology**

Hypercalcemia can be caused by hyperparathyroidism. Parathyroid hyperplasia, adenoma, or cancer can cause hyperparathyroidism. It is important for the surgeon planning resection of the parathyroid abnormality to know how many parathyroid glands are involved and their locations. Preoperative parathyroid scanning is the most accurate method of providing this information. Parathyroid hyperplasia causes enlargement of all four parathyroid glands. A parathyroid adenoma or cancer, however, causes enlargement of only one parathyroid gland and suppression of the other three glands.

Parathyroids are located most commonly on the lateral borders of the thyroid lobes—two on each side. However, parathyroid anatomic location varies considerably, and they may be located anywhere from the upper neck to the lower mediastinum. Parathyroid scanning is also done immediately before surgery to help the surgeon identify the parathyroid glands and particularly the pathologic glands. In this test, the scan is performed on the parathyroid glands as described previously. In the operating room, if the preoperative scan result is abnormal, the surgeon scans the suspect area of the neck with a hand-held gamma ray detector probe. Increased counts are noted in the regions where the parathyroid abnormalities are located.

**Contraindications**

- Patients who are pregnant

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Tell the patient that fasting is usually not required.

**During**

- Note the following procedural steps for the single tracer double phase (STDP) method.
  1. Technetium-99m sestamibi is injected intravenously.

## 678 parathyroid scan

2. At 15 minutes and 3 hours, the patient is placed in a supine position. For planar images, the detector is passed over the neck and upper chest area, and the radioactive counts are recorded and displayed.
  3. For SPECT/ computed tomography (CT) screening, the patient is placed in the appropriate unit in the supine position.
  4. Initially, the tracer lights up both the thyroid and parathyroid glands. At 3 hours, the tracer remains only in the pathologic parathyroid tissue.
- Note that this study is performed by a nuclear technologist and interpreted by a physician in the nuclear medicine department.
- PT** Tell the patient that no discomfort is associated with this study.

### After

- PT** Assure the patient that the dose of radioactive technetium used in this test is minute and therefore harmless.

### Abnormal findings

Parathyroid adenoma, carcinoma, or hyperplasia

Aberrantly placed parathyroid tissue in the upper neck, thyroid gland, or mediastinum

### notes

## Parkinson disease testing

**Type of test** Microscopic

### Normal findings

No peripheral synucleinopathy

### Test explanation and related physiology

Synucleinopathy refers to any degenerative disease of the central nervous system in which there is an excessive accumulation of alpha-synuclein (a brain protein) in the neurons. The synucleinopathies include Parkinson disease (PD), dementia with Lewy bodies, and multiple system atrophy. This test is used to diagnose PD. Peripheral synucleinopathy is also being used by researchers who strive to enroll patients in treatment trials early in the course of the disease.

In PD, alpha-synuclein proteins are found not only in the brain but throughout the body. Brain tissue examinations will most accurately identify these proteins. The submandibular gland, however, is the easiest to reach and potentially the safest, posing the fewest risks. The tissue is stained for phosphorylated alpha-synuclein. Rich staining indicates PD. Unfortunately, this test is only about 75% sensitive in early PD subjects. False positives may occur but may represent prodromal PD.

The DaT scan (Ioflupane I-123) is a radiopharmaceutical that can visualize striatal dopamine transporters as an adjunct to other diagnostic tests for PD. After intravenous (IV) injection, Ioflupane I-123 reversibly binds to the presynaptic dopamine transporter (DaT) in the striatum and is visualized using single photon emission computed tomography (SPECT). DaT is reduced 50% to 70% in patients with PD. When combined with olfactory testing such as the University of Pennsylvania Smell Identification Test (UPSIT-40) (because most patients with PD lose some sense of smell), the diagnosis of PD can be more strongly confirmed with fewer patients falsely diagnosed as having PD.

### Potential complications

- Submandibular gland bleeding, swelling, and inflammation

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Tell the patient no fasting is required

## 680 parkinson disease testing

### During

- Local anesthesia is used to numb the area.
- An otolaryngologist inserts a needle through the skin into the submandibular gland. Three to six needle core biopsies are obtained.

### After

- A dressing is applied to the neck area.
- Check the site for bleeding.

### Abnormal findings

PD

notes

**partial thromboplastin time, activated** (APTT, Partial thromboplastin time [PTT])**Type of test** Blood**Normal findings**

APTT: 30-40 seconds

PTT: 60-70 seconds

Patients receiving anticoagulant therapy: 1.5-2.5 times the control value in seconds

**Possible critical values**

APTT: &gt; 70 seconds

PTT: &gt; 100 seconds

**Test explanation and related physiology**

The PTT test is used to assess the intrinsic system and the common pathway of clot formation. PTT evaluates factors I (fibrinogen), II (prothrombin), V, VIII, IX, X, XI, and XII. When any of these factors exist in inadequate quantities, such as in hemophilia A and B or consumptive coagulopathy, PTT is prolonged. Because factors II, IX, and X are vitamin K-dependent factors, biliary obstruction, which precludes GI absorption of fat and fat-soluble vitamins (e.g., vitamin K), can reduce their concentration and thus prolong PTT. Because coagulation factors are made in the liver, hepatocellular diseases will also prolong PTT.

*Heparin* has been found to inactivate prothrombin (factor II) and prevent the formation of thromboplastin. These actions prolong the intrinsic clotting pathway for approximately 4 to 6 hours after each dose of heparin. Thus heparin is capable of providing therapeutic anticoagulation. The appropriate dose of heparin can be monitored by PTT. Test results are given in seconds, along with a control value. The control value may vary slightly from day to day because of the reagents used.

Recently activators have been added to PTT test reagents to shorten normal clotting time and provide a narrow normal range. This shortened time is called the *activated* PTT (APTT). Normal APTT is 30 to 40 seconds. Desired ranges for therapeutic anticoagulation are 1.5 to 2.5 times normal (e.g., 70 seconds). The APTT specimen should be drawn 30 to 60 minutes before the patient's next heparin dose is given. If the APTT is less than 50 seconds, the patient may not be receiving therapeutic

## 682 partial thromboplastin time, activated

anticoagulation and need more heparin. An APTT greater than 100 seconds indicates that too much heparin is being given; the risk of serious spontaneous bleeding exists when the APTT is this high. The effects of heparin can be reversed immediately by the administration of protamine sulfate. Anti-Xa (see p. 77) can also be used to monitor heparin dosing and determine dosing of protamine sulfate.

APTT is also used to determine *activated protein C (APC) resistance* (see p. 397). APC resistance testing is performed in the evaluation of thrombotic patients. A standard APTT test is performed first in the absence of, and then in the presence of, commercially available APC. The APTT is normally prolonged in the presence of APC because of APC's anticoagulant inhibiting coagulation factors VIII and V. An abnormality is detected if the APTT is not prolonged. This results from a resistance to APC, most commonly (80%) caused by an abnormal factor V Leiden.

Normal pregnancy has been associated with prolonged APTT times. Factor XI deficiency and antiphospholipid antibody are two major abnormalities identified in pregnant patients with prolonged APTT. These coagulopathies are not usually associated with excessive bleeding or thromboembolism. At the same time, pregnancy or oral contraceptive (OC) use is also associated with prolonged APTT. This is due, in large part, to a factor V mutation (factor V Leiden), which potentiates the prothrombotic effect of OC.

### Interfering factors

- ▼ Drugs that may prolong PTT test values include antihistamines, ascorbic acid, chlorpromazine, heparin, and salicylates.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: blue
- If the patient is receiving heparin by intermittent injection, plan to draw the blood specimen for the APTT 30 minutes to 1 hour before the next dose of heparin.
- If the patient is receiving continuous heparin, draw the blood at any time.
- Apply pressure to the venipuncture site. Remember that if the patient is receiving anticoagulants or has coagulopathies, the bleeding time will be increased.

- Assess the patient to detect possible bleeding. Check for blood in the urine and all other excretions and assess the patient for bruises, petechiae, and low back pain.
- If severe bleeding occurs, note that the anticoagulant effect of heparin can be reversed by parenteral administration of protamine sulfate.

### Abnormal findings

#### ▲ Increased levels

Acquired or congenital clotting factor deficiencies (e.g., hypofibrinogenemia, von Willebrand disease, and hemophilia)  
Cirrhosis of the liver  
Vitamin K deficiency  
Disseminated intravascular coagulation (DIC)  
Heparin administration

#### ▼ Decreased levels

Early stages of DIC  
Extensive cancer

### notes

## parvovirus B19 antibody

**Type of test** Blood

### Normal findings

Negative for immunoglobulin M (IgM)- and IgG-specific antibodies to parvovirus B19

### Test explanation and related physiology

The parvovirus B19 is known to be a human pathogen. Many of the severe manifestations of B19 viremia relate to the ability of the virus to infect and lyse RBC precursors in the bone marrow.

Erythema infectiosum is the most common manifestation of parvovirus B19 infection and occurs predominantly in children. This pathogen is also referred to as *fifth disease* or *academy rash*. Serologic testing for parvovirus B19-specific IgM and IgG antibodies can be done. Acute infections can be determined by parvovirus B19-compatible symptoms and the presence of IgM antibodies that remain detectable up to a few months. Past infection or immunity is documented by IgG antibodies that persist with IgM antibodies. Fetal infection may be recognized by hydrops fetalis and the presence of B19 DNA in amniotic fluid or fetal blood.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: verify with laboratory
- PT** Inform the patient that it normally requires approximately 2 to 3 days to get test results.

### Abnormal findings

#### ▲ Increased levels

- Erythema infectiosum (fifth disease)
- Joint arthralgia and arthritis
- Hydrops fetalis
- Fetal loss
- Transient aplastic anemia
- Chronic anemia in immunodeficient patients
- Bone marrow failure

notes



## pelvic ultrasonography (Pelvic ultrasonography in pregnancy, Obstetric ultrasonography, Vaginal ultrasound)

**Type of test** Ultrasound

### Normal findings

Normal fetal and placental size and position

### Test explanation and related physiology

Ultrasound examination of the female patient is a harmless, noninvasive method of evaluating the female genital tract and fetus. Pelvic ultrasonography can be performed with the transducer placed on the anterior abdomen (see [Figure 1](#), p. 1) or in the vagina with a specially designed vaginal probe. Vaginal ultrasound adds significant accuracy in identifying paracervical, endometrial, and ovarian pathology that otherwise may not be detected with the anterior abdominal probe.

Pelvic ultrasonography may be useful in the *obstetric patient* in the following circumstances:

- Making an early diagnosis of normal pregnancy and abnormal pregnancy (e.g., tubal pregnancy)
- Identifying multiple pregnancies
- Differentiating a tumor from a normal pregnancy
- Determining the age of the fetus by the diameter of the head
- Measuring the rate of fetal growth
- Identifying placental abnormalities
- Diagnosing ectopic pregnancy
- Providing a realistic image of the fetus using three- or four-dimensional imaging
- Evaluating the kidneys and upper collecting system
- Localizing the placenta before amniocentesis
- Making a differential diagnosis of various uterine and ovarian enlargements (e.g., polyhydramnios)

Ultrasonography is a very accurate and easily performed screening test to recognize risks of fetal abnormalities (see amniotic fluid index, p. 406). *Fetal nuchal translucency (FNT)* is an ultrasound measurement of subcutaneous edema in the neck region of the fetus. It is performed at 10 to 14 weeks of gestation. Major heart defects, trisomy 21, and other genetic defects are associated with increased edema in this location at this age of gestation. With FNT, these abnormalities can be identified earlier in the pregnancy, when abortion is still possible.

## 686 pelvic ultrasonography

Pelvic ultrasonography is used in the *nonpregnant woman* to monitor the endometrium in patients who take tamoxifen and to aid in the diagnosis of the following:

- Ovarian cyst or tumor
- Tubo-ovarian abscess
- Uterine fibroids or cancer
- Pelvic inflammatory disease (PID)
- Thickened uterine endometrium (stripe) (caused by cancer, hyperplasia, and so on)

Pelvic ultrasonography is widely used in the evaluation of infertility of the female. *Saline infusion sonography (SIS)* is a simple and inexpensive method to determine patency of the uterus and tubes. In SIS, saline is injected under ultrasound guidance. This provides better visualization of any intrauterine abnormalities. Air bubbles are then injected and imaged through the tubes to demonstrate patency.

When a woman is unable to visualize or palpate the string of an intrauterine device (IUD), ultrasound is indicated to determine whether the IUD has perforated the uterus, been expelled, or been incorporated with an intrauterine pregnancy.

### Contraindications

- Patients with latex allergy. Vaginal ultrasound requires placement of the probe in a latex condom–like sac.

### Interfering factors

- Patients who have had recent GI contrast studies because barium creates severe distortion of reflective sound waves
- Patients with air-filled bowels because gas does not transmit the sound waves well
- Failure to fill the bladder may make the image uninterpretable.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting or sedation is required.
- PT** Assure the patient that this study has no known deleterious effects on maternal or fetal tissues.
- PT** Give the patient 3 to 4 glasses (200-350 mL) of water or another liquid 1 hour before abdominal ultrasonography, and instruct the patient *not* to void until after the procedure is completed. This will permit visualization of the bladder, which is used as a reference point in pelvic anatomy.

- If a transabdominal ultrasonography is required urgently and there is not time to fill the bladder by ingestion of fluids, a bladder catheter is inserted and the bladder is filled with water.
- No water is required for vaginal ultrasonography.

### During

- Note the following procedural steps:
    1. In the ultrasound room, the patient is placed in the supine position on the examining table.
    2. The ultrasonographer applies a greasy conductive paste to the abdomen to enhance sound transmission and reception.
    3. A transducer is passed over the skin.
    4. If a vaginal probe is used, it is inserted via the vagina and angled to identify the various parts of the pelvis.
    5. Pictures are taken of the reflections.
    6. During the examination, fetal structures are usually pointed out to the mother.
  - Note that this procedure is performed in approximately 20 minutes.
- PT** Inform the patient that no discomfort is associated with this study other than having a full bladder and the urge to void.

### After

- Remove the lubricant from the patient's skin.

### Abnormal findings

Tubal pregnancy

Abdominal pregnancy

Hydatidiform mole of the uterus

Intrauterine growth retardation

Multiple pregnancy

Fetal death

Abruptio placentae

Abnormal fetal position

Fetal anomalies

Placenta previa

Polyhydramnios

Neoplasm of the ovaries, uterus, or fallopian tubes

Cysts

Abscesses

Hydrocephalus of the fetus

### notes

## pepsinogen

**Type of test** Blood/urine

### Normal findings

Pepsinogen I: 28-100 ng/ml

Pepsinogen II: < 22 ng/mL

### Test explanation and related physiology

Pepsinogens are secreted in the stomach and are made in the oxyntic gland mucosa of the proximal stomach. When exposed to gastric acid, pepsinogen is converted to pepsin, an active enzyme that is proteolytic and promotes digestion. Patients with gastric atrophy or pernicious anemia (PA) or those who have had a gastrectomy have low levels of pepsinogen I. Pepsinogen I levels are elevated in patients with ulcer disease. Pepsinogen I has been used as a subclinical marker of increased risk for stomach cancer. Pepsinogen I can also be measured in the urine.

Pepsinogen II is made by oxyntic gland mucosa cells that are in the distal stomach and proximal duodenum. Because PA generally affects the proximal stomach, diminished levels of pepsinogen I with normal levels of pepsinogen II are strongly supportive of PA.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red
- PT** Tell the patient that antacids or other medications affecting stomach acidity or GI motility should be discontinued, if possible, for at least 48 hours before collection. Verify with the laboratory or health-care provider.
- The specimen is usually frozen because most testing is performed in a central laboratory.

### Abnormal findings

#### ▼ Decreased values

- Pernicious anemia
- Gastric atrophy
- Chronic gastritis
- Peptic ulcer disease

#### notes

## pericardiocentesis

**Type of test** Fluid analysis

### Normal findings

Less than 50 mL of clear, straw-colored fluid without evidence of any bacteria, blood, or malignant cells

### Test explanation and related physiology

Pericardiocentesis, which involves the aspiration of fluid from the pericardial sac with a needle, may be performed for therapeutic and diagnostic purposes. Therapeutically, the test is performed to relieve cardiac tamponade by removing fluid and improving diastolic filling. Diagnostically, pericardiocentesis is performed to remove a sample of pericardial fluid for laboratory examination to determine the cause of the fluid. This is similar to the evaluation described for peritoneal and pleural fluid (see p. 671 and p. 864, respectively).

### Contraindications

- Patients who are uncooperative because of the risk of lacerations to the epicardium or coronary artery
- Patients with a bleeding disorder. Inadvertent puncture of the myocardium may create uncontrollable bleeding, leading to tamponade.

### Potential complications

- Laceration of the coronary artery or myocardium
- Needle-induced ventricular arrhythmias
- Myocardial infarction
- Pneumothorax caused by inadvertent puncture of the lung
- Liver laceration caused by inadvertent puncture
- Pleural infection
- Vasovagal arrest

### Procedure and patient care

#### Before


**PT** Explain the procedure to the patient.

- Obtain informed consent for this procedure.
- Restrict fluid and food intake for 4 hours (if possible).
- Obtain IV access for infusion of fluids and cardiac medications if required.
- Administer pretest medication. Atropine may be given to prevent the vasovagal reflex of bradycardia and hypotension.

## 690 pericardiocentesis

### During

- Note the following procedural steps:
  1. The patient is placed in the supine position.
  2. An area in the fifth to sixth intercostal space at the left sternal margin (or subxiphoid) is prepared and draped.
  3. After skin anesthesia is performed, a large-bore pericardiocentesis needle is placed on a 50-mL syringe and introduced into the pericardial sac ([Figure 33](#)).
  4. An electrocardiographic lead is often attached by a clip to the needle to identify any ST-segment elevations.
  5. Pericardial fluid is aspirated.
  6. Some patients who have recurring cardiac tamponade may require placement of an indwelling pericardial catheter.

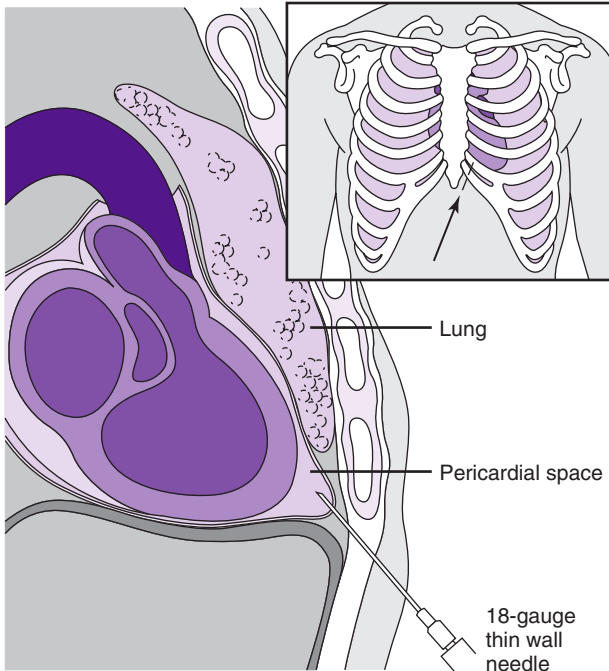

**FIG. 33** Pericardiocentesis using subxiphoid route for aspiration of pericardial fluid. The 18-gauge needle is introduced at a 30- to 40-degree angle.

7. With certain types of pericarditis, medications (e.g., antibiotics, antineoplastic drugs, corticosteroids) may be instilled during pericardiocentesis.
- Note that a physician usually performs this procedure in the cardiac catheterization laboratory, operating room, or emergency department in approximately 10 to 20 minutes.
- PT** Tell the patient that this procedure is associated with very little discomfort.

### After

- Closely monitor the patient's vital signs. An increased temperature may indicate infection. Pericardial bleeding would be marked by hypotension or pulsus paradoxus.
- Check the dressing frequently for drainage.
- Label and number the specimen tubes that contain the pericardial fluid and deliver them to the appropriate laboratories for examination. Note the following possibilities:
  1. Usually the fluid is taken to the chemistry laboratory, where the color, turbidity, glucose, albumin, protein, and lactic dehydrogenase levels are obtained.
  2. A tube of blood often goes to the hematology laboratory, where red and white blood cells are evaluated.
  3. The bacteriology laboratory performs routine cultures, Gram stains, fungal studies, acid-fast bacilli smears, and so on.
  4. When malignancy is suspected, the fluid should be sent for cytology.
- Apply a sterile dressing to the catheter if one has been left for continuing pericardial drainage.
- Establish a closed system if continued pericardial drainage is required.
- Note that to minimize infection, pericardial catheters, if used, are usually removed after 2 days.

### Abnormal findings

Pericarditis	Metastatic cancer
Uremia	Blunt or penetrating cardiac trauma
Hypoproteinemia	Rupture of ventricular aneurysm
Congestive heart failure	Collagen vascular diseases

### notes

**pheochromocytoma suppression and provocative testing** (Clonidine suppression test [CST], Glucagon stimulation test)**Type of test** Blood**Normal findings**

Glucagon stimulation

Norepinephrine: &lt; 3 times basal levels

Clonidine suppression

Norepinephrine: &gt; 50% reduction in basal levels or &lt; 500 pg/mL

Epinephrine: &gt; 50% reduction in basal levels or &lt; 275 pg/mL

**Test explanation and related physiology**

In patients with significantly high blood pressure refractory to treatment, the diagnosis of pheochromocytoma is often considered. When catecholamine levels are excessive (norepinephrine > 2000 pg/mL), the diagnosis is easily made. However, when basal levels are not significantly elevated, it is difficult to differentiate essential hypertension from a functioning pheochromocytoma. Suppression and provocative tests may be necessary. Glucagon is used as the provocative agent. In patients with pheochromocytoma, the response is accentuated. Clonidine is normally a potent suppressor of catecholamine production, yet it has little to no effect on catecholamines in patients with pheochromocytoma. Testing of metanephrines (see p. 621) provides higher diagnostic sensitivity than catecholamine assays in screening for pheochromocytoma.

**Contraindications**

- Hypovolemic/dehydrated patients. These patients should not have suppression testing because they could experience a precipitous drop in blood pressure.

**Potential complications**

- Drowsiness during CST
- Hypotension during CST, especially in patients treated aggressively for hypertension
- Extremely high blood pressure during provocative testing

**Interfering factors**

- False suppression with CST may occur in patients with low basal catecholamine levels.



## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- Identify the medications being administered before testing.
  - The patient must be reclining calmly for 30 minutes before testing.

### During

- Collect a venous blood sample from an antecubital vein in a heparinized tube for determination of basal catecholamine levels.
- Monitor vital signs closely throughout the testing period.

#### *Glucagon provocative test*

- Administer a prescribed dose of glucagon intravenously.
- Two minutes later, obtain a blood specimen as described above.

#### *Clonidine suppression test*

- Administer a prescribed dose of clonidine orally.
- Three hours later, obtain a blood specimen as described earlier.

### After

- Monitor vital signs for at least 1 hour after conclusion of the procedure.

## Abnormal findings

Pheochromocytoma

notes

**phosphate** ( $\text{PO}_4$ , Phosphorus [P])**Type of test** Blood**Normal findings**

Adult: 3.0-4.5 mg/dL or 0.97-1.45 mmol/L (SI units)

Elderly: values slightly lower than adult

Child: 4.5-6.5 mg/dL or 1.45-2.1 mmol/L (SI units)

Newborn: 4.3-9.3 mg/dL or 1.4-3 mmol/L (SI units)

**Possible critical values**

&lt; 1 mg/dL

**Test explanation and related physiology**

Phosphorus in the body exists in the form of a phosphate. The terms *phosphorus* and *phosphate* are used interchangeably throughout this and other discussions. Most of the phosphate in the body is a part of organic compounds. Only a small part of total body phosphate is inorganic phosphate (i.e., not part of another organic compound). It is the *inorganic* phosphate that is measured when one requests a phosphate, phosphorus, inorganic phosphorus, or inorganic phosphate. Most of the body's inorganic phosphorus is intracellular and combined with calcium within the skeleton; however, approximately 15% of the phosphorus exists in the blood as a phosphate salt.

Dietary phosphorus is absorbed in the small bowel. The absorption is very efficient, and only rarely is hypophosphatemia caused by GI malabsorption. Phosphorus levels are determined by calcium metabolism, parathormone (PTH), renal excretion, and, to a lesser degree, intestinal absorption. Because an inverse relationship exists between calcium and phosphorus, a decrease of one mineral results in an increase in the other. The regulation of phosphate by PTH is such that PTH tends to decrease phosphate resorption in the kidney. PTH and vitamin D, however, tend to stimulate phosphate absorption weakly within the gut. Hypophosphatemia may have four general causes: shift of phosphate from extracellular to intracellular, renal phosphate wasting, loss from the GI tract, and loss from intracellular stores. Hyperphosphatemia is usually secondary to increased intake or an inability of the kidneys to excrete phosphate.

## Interfering factors

- Laxatives or enemas containing sodium phosphate can increase phosphorus levels.
- Recent carbohydrate ingestion, including IV glucose administration, causes decreased phosphorus levels because phosphorus enters the cell with glucose.
- ▼ Drugs that may cause *increased* levels include methicillin and vitamin D (excessive).
- ▼ Drugs that may cause *decreased* levels include albuterol, anesthesia agents, antacids, estrogens, insulin, mannitol, and oral contraceptives.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes
- Blood tube commonly used: red
- If indicated, discontinue IV fluids with glucose for several hours before the test.
- Avoid hemolysis. Handle the tube carefully.
- Use a heel stick to draw blood from infants.

## Abnormal findings

▲ **Increased levels (hyperphosphatemia)**  
 Renal failure  
 Increased dietary or IV intake of phosphorus  
 Acromegaly  
 Hypoparathyroidism  
 Bone metastasis  
 Sarcoidosis  
 Hypocalcemia  
 Liver disease  
 Acidosis  
 Rhabdomyolysis  
 Advanced lymphoma or myeloma  
 Hemolytic anemia

▼ **Decreased levels (hypophosphatemia)**  
 Inadequate dietary ingestion of phosphorus  
 Chronic antacid ingestion  
 Hyperparathyroidism  
 Hypercalcemia  
 Chronic alcoholism  
 Vitamin D deficiency  
 Diabetic acidosis  
 Hyperinsulinism  
 Rickets (childhood)  
 Osteomalacia (adult)  
 Malnutrition  
 Alkalosis  
 Sepsis (gram negative)

notes

**phospholipase A2 receptor antibodies (anti-pla2r)****Type of test** Blood**Normal findings**

Negative: &lt; 14 RU/mL

Borderline: 14-19 RU/mL

Positive:  $\geq$  20 RU/mL**Test explanation and related physiology**

Membranous nephropathy (MN) is a rare disease in which immune complexes deposit at the glomerular basement membrane, causing damage to the filtration barrier, resulting in proteinuria and nephrotic syndrome. Nearly 70% of patients with primary MN (pMN) have autoantibodies directed against the phospholipase A2 receptor (PLA2R). There is also evidence that levels of anti-PLA2R autoantibodies correlate well with disease activity and progression. The presence of anti-PLA2R antibodies can be used to differentiate pMN from other causes of nephrotic syndrome (or secondary membranous nephropathy). Among patients awaiting kidney transplantation, higher levels of anti-PLA2R could predict those more likely for recurrence of pMN after transplantation. Furthermore, successful treatment of pMN will result in reduced PLA2R antibodies. An increase in antibodies would indicate a worsening of the disease.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

**Abnormal findings**

Primary membranous nephropathy (pMN)

**notes**

## PI-linked antigen (Phosphatidylinositol antigen)

**Type of test** Blood

### Normal findings

RBCs

Type I (normal expression): 99%-100%

Type II (partial deficient): 0.00%-0.99%

Type III (deficient): 0.00%-0.01%

Granulocytes: 0.00%-0.01%

Monocytes: 0.00%-0.05%

### Test explanation and related physiology

*Paroxysmal nocturnal hemoglobinuria (PNH)* is an acquired hematologic disorder characterized by nocturnal hemoglobinuria, chronic hemolytic anemia, thrombosis, pancytopenia, and, in some patients, acute or chronic myeloid malignancies. PNH appears to be a hematopoietic stem cell disorder that affects erythroid, granulocytic, and megakaryocytic cell lines. The abnormal cells in PNH have been shown to lack glycosylphosphatidylinositol (GPI)-linked proteins in RBCs and WBCs. Mutations in the *phosphatidylinositol glycan A (PIGA)* gene have been identified consistently in patients with PNH, thus confirming the biological defect in this disorder.

Flow cytometric immunophenotyping of peripheral blood (WBC and RBC) is performed to evaluate the presence or absence of PI-linked antigens (CD14, FLAER, and/or CD59 antigens), using monoclonal antibodies directed against them. These proteins are absent on the blood cells of patients with PNH. Determination of PI-linked antigens is not only useful in the diagnosis of PNH but also in monitoring of the disease.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: yellow

### Abnormal findings

#### ▼ Decreased levels

PNH

notes

## placental growth factor (PGF)

**Type of test** Blood

### Normal findings

Nonpregnant women: < 50 pg/mL

Pregnant women at 22 weeks' gestation:  $\geq$  200 pg/mL

### Test explanation and related physiology

Preeclampsia is one of the most common medical complications of pregnancy and is associated with considerable maternal and neonatal morbidity and mortality. PGF is a protein that is produced during pregnancy by the placental trophoblast. Normally, PGF rises steadily throughout pregnancy to levels exceeding 500 pg/mL. When PGF is tested at 13 to 16 weeks of gestation and it is found to be significantly decreased, the patient is at considerable risk for preeclampsia. Likewise, patients with markedly increased levels of soluble fms-like tyrosine kinase-1 (sFlt-1), which is a known inhibitor of PGF, are also at risk for preeclampsia.

Preeclampsia can be predicted by a combination of factors in the maternal history, including black racial origin, high body mass index, family history of preeclampsia, and personal history of preeclampsia. Screening at risk women with PGF, sFlt-1, PAPP-A (see p. 730), and uterine artery Doppler would identify women who would develop early and late preeclampsia. Intensive maternal and fetal monitoring could improve pregnancy outcomes.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▼ Decreased blood levels

Hypertension  
Preeclampsia

notes

**plasminogen** (Fibrinolysin)**Type of test** Blood**Normal findings**

2.4-4.4 Committee on Thrombolytic Agents (CTA) units/mL

**Test explanation and related physiology**

This test is used to diagnose suspected plasminogen deficiency in patients who have multiple thromboembolic episodes. Plasminogen is a protein involved in the fibrinolytic process of intravascular blood clot dissolution (see [Figure 10](#), p. 251). Plasminogen is converted to plasmin by proteolytic cleavage. Plasmin can destroy fibrin and dissolve clots.

Plasminogen levels are occasionally measured during fibrinolytic therapy (for coronary and peripheral arterial occlusion) and are diminished with full fibrinolysis. Decreased levels of plasminogen are also found in hyperfibrinolytic states. Because plasminogen is made in the liver, patients with cirrhosis or other severe liver diseases can be expected to have decreased levels. Inflammatory conditions may have mild elevations of plasminogens.

**Interfering factors**

- Pregnancy and eclampsia are associated with elevations.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: blue

**Abnormal findings****▲ Increased levels**

Pregnancy  
Inflammatory conditions

**▼ Decreased levels**

Hyperfibrinolytic state  
(e.g., DIC, fibrinolysis)  
Primary liver disease  
Syndrome associated with hypercoagulation  
Rare congenital deficiencies  
Malnutrition

notes

**plasminogen activator inhibitor 1 (PAI-1)****Type of test** Blood**Normal findings**

Antigen assay: 2-46 ng/mL

Activity: &lt; 31.1 IU/mL

**Test explanation and related physiology**

PAI-1 is a protein that inhibits plasminogen activators. During fibrinolysis, tissue plasminogen activator (tPA) converts plasminogen into plasmin. Plasmin plays a critical role in fibrinolysis by degrading fibrin (see [Figure 10](#), p. 251). PAI-1 is the primary inhibitor of tPA and urokinase plasminogen activator (uPA) in the blood. PAI-1 limits the production of plasmin and keeps fibrinolysis in check.

Elevated levels of PAI-1 are associated with a predisposition to thrombosis, including veno-occlusive disease after bone marrow transplantation or high-dose chemotherapy. Familial thrombosis has been associated with inherited elevation of plasma PAI-1 activity. Increased levels of PAI-1 have also been reported in a number of conditions, including malignancy, liver disease, the postoperative period, septic shock, the second and third trimesters of pregnancy, obesity, coronary heart disease, and restenosis after coronary angioplasty. Increased levels may reduce the effectiveness of antithrombotic therapy. Patients with insulin resistance syndrome and diabetes mellitus tend to have increased PAI-1 levels.

Low plasma levels of the active form of PAI-1 have been associated with abnormal, clinically significant bleeding. Complete deficiency of PAI-1, either congenital or acquired, is associated with bleeding manifestations that include hemarthroses, hematomas, menorrhagia, easy bruising, and postoperative hemorrhage.

**Interfering factors**

- Because PAI-1 is an acute-phase reactant, it can become transiently elevated by infection, inflammation, or trauma.
- Levels increase during pregnancy.
- PAI-1 has a circadian rhythm, with the highest concentrations occurring in the morning and the lowest concentrations in the afternoon and evening.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no



- Blood tube commonly used: light blue
- Discard the first several milliliters of blood if PAI-1 is the only test being drawn. If multiple tests are being drawn, fill the light blue-top tube after any red-top tube.
- Gently invert the blood tube several times after collection.

### Abnormal findings

#### ▲ Increased levels

Acute coronary syndrome  
Coronary artery disease  
Restenosis after coronary  
angioplasty  
Infection  
Inflammation  
Trauma  
Diabetes mellitus  
Insulin resistance syndrome  
Pregnancy

#### ▼ Decreased levels

Bleeding disorders

notes

## platelet aggregation test (Ristocetin test)

**Type of test** Blood

### Normal findings

Dependent on the platelet agonist used

### Test explanation and related physiology

Platelet aggregation is an important part of hemostasis. Surrounding an area of acute blood vessel endothelial injury is a clump of platelets. Normal platelets adhere to this area of injury, and through a series of chemical reactions, they attract other platelets to the area. This is platelet aggregation, the first step of hemostasis. After this step, the normal coagulation factor water-fall occurs (see [Figure 10](#), p. 251).

Certain diseases that affect either platelet number or function can inhibit platelet aggregation and thereby prolong bleeding times. Congenital syndromes, uremia, myeloproliferative disorders, and certain drugs are also associated with abnormal platelet aggregation. If blood is passed through a heart–lung or dialysis pump, platelet injury can occur and aggregation can be reduced.

Ristocetin is most commonly used to induce platelet agglutination. The agglutination of the patient's blood with the addition of this type of product can help differentiate the disease that may be affecting platelet aggregation.

### Interfering factors

- Factors that may cause *increased* platelet aggregation include blood storage temperature, hyperbilirubinemia, hemoglobinemia, hyperlipidemia, and platelet count.
- ✎ Drugs that may cause *decreased* platelet aggregation include antibiotics, aspirin, nonsteroidal antiinflammatory agents, and thienopyridine antiplatelet drugs (e.g., ticlopidine, clopidogrel).

### Procedure and patient care

- See [inside front cover for Routine Blood Testing](#).
- Fasting: no
- Blood tube commonly used: blue
- Remember that abnormalities in platelet aggregation can prolong bleeding time, and a significant hematoma at the venipuncture site may occur.

## Abnormal findings

Hypoactive platelet aggregation

Various congenital disorders (e.g., Wiskott–Aldrich syndrome, Bernard–Soulier syndrome, glycogen storage, some types of von Willebrand disease)

notes

**platelet antibody detection** (Antiplatelet antibody detection)**Type of test** Blood**Normal findings**

No antiplatelet antibodies identified

**Test explanation and related physiology**

Immune-mediated destruction of platelets may be caused either by autoantibodies directed against antigens located on the same person's platelets or by alloantibodies that develop after exposure to transfused platelets received from a donor. These antibodies are usually directed to an antigen on the platelet membrane, such as *human lymphocyte antigen (HLA)* (see p. 518) or *platelet-specific antigens*, such as *PLA1* and *PLA2*.

Antibodies directed to platelets cause early destruction of the platelets and subsequent thrombocytopenia. Immunologic thrombocytopenia includes the following:

- *Idiopathic thrombocytopenia purpura (ITP)*. Platelet-associated IgG antibodies are detected in 90% of these patients.
- *Posttransfusion purpura*. This is usually associated with an antibody to ABO, HLA, or PLA antigens on the RBC.
- *Maternal-fetal platelet antigen incompatibility (neonatal thrombocytopenia)*. This occurs when the fetal platelet contains a PLA1 antigen that is absent in the mother. Neonatal thrombocytopenia can also occur if the mother has ITP autoantibodies that are passed through the placenta and destroy the fetal platelets.
- *Drug-induced thrombocytopenia*. Although a host of drugs are known to induce autoimmune-mediated thrombocytopenia, heparin is the most common and causes heparin-induced thrombocytopenia (HIT). There are two types of HIT, type I and type II, that may develop. Type I HIT is generally considered a benign condition and is not antibody-mediated. In type II HIT, thrombocytopenia is usually more severe and is antibody-mediated. Type II HIT is caused by an IgG antibody and usually occurs after 6 to 8 days of IV heparin therapy. The diagnosis of HIT is confirmed by identifying *heparin-induced thrombocytopenia antibodies (HITAs)*.

Other drugs known to cause antiplatelet antibodies include analgesics (salicylates), antibiotics (cephalosporins), cimetidine, diuretics, heavy metals (e.g., gold), hypnotics, oral hypoglycemic agents, quinidine-like drugs, and many others (e.g., digoxin).

## Interfering factors

- Blood transfusion may cause the development of isoantibodies to HLA antigens on the platelets or RBCs.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

## Abnormal findings

### ▲ Increased levels

Idiopathic thrombocytopenia purpura

Neonatal thrombocytopenia

Posttransfusion purpura

HIT

Drug-induced thrombocytopenia

Paroxysmal hemoglobinuria

notes

## platelet count (Thrombocyte count)

**Type of test** Blood

### Normal findings

Adult/elderly: 150,000-400,000/mm<sup>3</sup> or 150-400 × 10<sup>9</sup>/L  
(SI units)

Premature infant: 100,000-300,000/mm<sup>3</sup>

Newborn: 150,000-300,000/mm<sup>3</sup>

Infant: 200,000-475,000/mm<sup>3</sup>

Child: 150,000-400,000/mm<sup>3</sup>

### Possible critical values

< 50,000 or > 1 million/mm<sup>3</sup>

### Test explanation and related physiology

The platelet count is a count of the number of platelets (thrombocytes) per cubic milliliter of blood. It is performed on patients who develop petechiae (small hemorrhages in the skin), spontaneous bleeding, or increasingly heavy menses. It is also used to monitor the course of the disease or therapy for thrombocytopenia or bone marrow failure.

Platelet activity is essential to blood clotting. Counts of 150,000 to 400,000/mm<sup>3</sup> are considered normal. Counts of less than 100,000/mm<sup>3</sup> are generally considered to indicate *thrombocytopenia*; *thrombocytosis* (*thrombocythemia*) is generally said to exist when counts are greater than 400,000/mm<sup>3</sup>. Common associative diseases with thrombocytosis are iron-deficiency anemia and malignancy (leukemia, lymphoma, or solid tumors, such as those involving the colon). Thrombocytosis may also occur with polycythemia vera, postsplenectomy syndromes, and a variety of acute/chronic infections or inflammatory processes.

Causes of thrombocytopenia (decreased number of platelets) include the following:

- Reduced production of platelets (secondary to bone marrow failure or infiltration of fibrosis, tumor, and so on)
- Sequestration of platelets (secondary to hypersplenism)
- Accelerated destruction of platelets (secondary to antibodies [see antiplatelet antibodies, p. 704], infections, drugs, prosthetic heart valves)
- Consumption of platelets (secondary to DIC)

- Platelet loss from hemorrhage
- Dilution with large volumes of blood transfusions that contain very few, if any, platelets

### Interfering factors

- Living at high altitudes may cause increased platelet levels.
- Because platelets can clump together, automated counting is subject to at least a 10% to 15% error.
- Strenuous exercise may cause increased levels.
- Decreased levels may be seen before menstruation.
- ✚ Drugs that may cause *increased* levels include oral contraceptives.
- ✚ Drugs that may cause *decreased* levels include chemotherapeutic agents, chloramphenicol, colchicine, H<sub>2</sub>-blocking agents (cimetidine, ranitidine), hydralazine, indomethacin, isoniazid, levofloxacin, quinidine, streptomycin, sulfonamides, thiazide diuretics, and tolbutamide.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender
- If the results indicate that the patient has a serious platelet deficiency, perform the following steps:
  1. Observe the patient for signs and symptoms of bleeding.
  2. Check for blood in the urine and all excretions.
  3. Assess the patient for bruises, petechiae, bleeding from the gums, epistaxis, and low back pain.
  4. Reassess all venipuncture sites for signs of hematoma formation.

### Abnormal findings

▲ **Increased levels  
(thrombocytosis)**

- Malignant disorder
- Polycythemia vera
- Postsplenectomy  
  syndrome
- Rheumatoid arthritis
- Iron-deficiency anemia

▼ **Decreased levels  
(thrombocytopenia)**

- Hypersplenism
- Hemorrhage
- Immune  
  thrombocytopenia
- Leukemia and other  
  myelofibrosis disorders
- Thrombotic  
  thrombocytopenia
- Inherited  
  thrombocytopenia  
  disorders (e.g.,  
  Wiskott–Aldrich,  
  Bernard–Soulier,  
  or Zieve syndromes)
- DIC
- Systemic lupus  
  erythematosus
- Pernicious anemia
- Some hemolytic anemias
- Cancer chemotherapy
- Acute or chronic infection

notes



## platelet function assay (Platelet closure time, PCT, Aspirin resistance tests, Bleeding time [BT], 11-Dehydro-thromboxane B2)

**Type of test** Blood, urine

### Normal findings

Platelet closure time (blood)

Collagen/adenosine-5'-diphosphate (CADP): 64-120 seconds

Collagen/epinephrine (CEPI): 89-193 seconds

11-Dehydro-thromboxane B2 (urine)

Males: 0-1089 pg/mg of creatinine

Females: 0-1811 pg/mg of creatinine

Bleeding time (blood)

1-9 minutes (Ivy method)

### Test explanation and related physiology

Platelet dysfunction may be acquired, inherited, or induced by platelet-inhibiting agents. It is clinically important to assess platelet function as a potential cause of a bleeding diathesis (epistaxis, menorrhagia, postoperative bleeding, or easy bruising). The most common causes of platelet dysfunction are related to uremia, liver disease, von Willebrand disease (vWD), and exposure to such agents as acetyl salicylic acid (ASA, aspirin). Several tests are used to evaluate platelet function. To measure *bleeding time*, a small standard superficial incision is made in the forearm, and the time required for the bleeding to stop is recorded. If a larger skin vessel is lacerated during the test, the bleeding time will be artificially prolonged. A repeat test is required. Compared with other alternatives, bleeding time (BT) is a bit more labor intensive, and its accuracy is heavily dependent on operator skills. Furthermore, its results are not easily reproduced and quantified. The platelet aggregation study (see p. 702) may also have similar problems. With the development of an automated platelet function analyzer device, clinical laboratories can easily measure *platelet closure time* (PCT) to quantify platelet function. Furthermore, PCT can differentiate aspirin effects from other causes of platelet dysfunction.

This test can also be used to determine resistance of aspirin's therapeutic anticoagulation effects on platelets. This is one of several *aspirin resistance tests* that are performed to determine the effectiveness of aspirin on inhibiting platelet aggregation and thereby protecting the patient from vascular thromboembolic disease (Table 31).

TABLE 31 Platelet closure time

	Normal	ASA effect	Intrinsic platelet disorders
CEPI membrane	Normal	Abnormal	Abnormal
CADP membrane	Normal	Normal	Abnormal

ASA, Aspirin; CADP, collagen/adenosine-5'-diphosphate; CEPI, collagen/epinephrine.

Aspirin resistance can be determined by platelet closure time or by measurement of *11-dehydro-thromboxane B2 (11-dTXB2)* in the urine. Thromboxane A2 is produced by the enzyme cyclo-oxygenase-1 (COX1) by activated platelets and still further stimulates platelet activation, platelet aggregation, and vasoconstriction. 11-dTXB2 is the stable, inactive metabolite of thromboxane A2. Urinary 11-dTXB2, therefore is an indication of platelet activation and aggregation. Elevated values are associated with an increased risk of acute ischemic stroke and myocardial infarction. Effective aspirin therapy should reduce the level of this metabolite in the urine. If not, the patient may be aspirin resistant and may be more safely treated with an alternative therapy such as increasing the dosage of aspirin or placing the patient on another antiplatelet medication.

Urinary 11-dTXB2 offers an advantage over blood aspirin resistance tests because it is not subject to interference from *in vitro* platelet activation caused by local vein trauma or insufficient anticoagulation during blood sample collection.

### Interfering factors

- Low hematocrit or platelet count can prolong PCT and BT.
- ✚ Aspirin and nonsteroidal antiinflammatory drugs (NSAIDs) can *prolong* PCT and DT test results or decrease urinary 11-dTXB2 levels.
- ✚ Thienopyridines can *prolong* test PCT.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: light blue
- Obtain a drug history to determine whether the patient has recently had aspirin or any other medications that may affect test results. Patient should not have taken NSAIDs within

72 hours or aspirin within 2 weeks before collection of a specimen for baseline analysis.

- For urinary 11-dTXB<sub>2</sub>, randomly collect 10 mL of urine. No preservative is necessary.

## Abnormal findings

### ▲ Prolonged times or increased values

#### *Intrinsic platelet defects*

Some myelodysplastic syndromes

Some myeloid leukemias

Some myeloproliferative neoplasms

Bernard–Soulier syndrome

Glanzmann thromboasthenia

Hermansky–Pudlak syndrome

Hereditary telangiectasia

#### *Platelet–blood vessel interaction defects*

von Willebrand disease

Cushing syndrome

Henoch–Schönlein syndrome

Uremia

#### *Elevated B<sub>2</sub> (11-dTXB<sub>2</sub>) (on aspirin therapy)*

Increased risk of thromboembolic disease

## notes

**platelet volume, mean** (Mean platelet volume [MPV])**Type of test** Blood**Normal findings**

7.4-10.4 fL

**Test explanation and related physiology**

The MPV is a measure of the volume of a large number of platelets determined by an automated analyzer. MPV is to platelets as mean corpuscular volume (see p. 772) is to the RBC.

The MPV varies with total platelet production. In cases of thrombocytopenia, despite normal reactive bone marrow (e.g., hypersplenism), the normal bone marrow releases immature platelets to attempt to maintain a normal platelet count. These immature platelets are larger, and the MPV is elevated. When bone marrow production of platelets is inadequate, the platelets that are released are small. This is reflected as a low MPV; in this way, the MPV is useful in the differential diagnosis of thrombocytopenic disorders.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender
- If the patient is known to have a low platelet count, perform the following steps:
  1. Observe the patient for signs and symptoms of bleeding.
  2. Check for blood in the urine and all excretions.
  3. Assess the patient for bruises, petechiae, bleeding of the gums, epistaxis, and low back pain.

**Abnormal findings****▲ Increased levels**

Valvular heart disease  
 Immune thrombocytopenia  
 Massive hemorrhage  
 B<sub>12</sub> or folate deficiency  
 Myelogenous leukemia

**▼ Decreased levels**

Aplastic anemia  
 Chemotherapy-induced  
 myelosuppression  
 Wiskott–Aldrich syndrome

**notes**

**plethysmography, arterial** (Ankle-brachial index [ABI])

**Type of test** Manometric

**Normal findings**

< 20 mm Hg difference in systolic blood pressure of the lower extremity compared with the upper extremity

Normal pulse wave amplitude showing a steep upswing; an acute, narrow peak; and a gentler downslope containing a dicrotic notch (normal arterial pulse wave)

Ankle/brachial ratio: 0.9-1.3

**Test explanation and related physiology**

Plethysmography is usually performed to rule out occlusive disease of the lower extremities; however, it also can identify arteriosclerotic disease in the upper extremities. This test requires one normal extremity against which the other extremities may be compared.

Arterial plethysmography is performed by applying three blood pressure cuffs to the proximal, middle, and distal parts of an extremity. Pressure readings are also taken in the upper arm (brachial) artery. These are then attached to a pulse volume recorder (plethysmograph), and each pulse wave can be displayed. A reduction in amplitude of a pulse wave in any of the three cuffs indicates arterial occlusion immediately proximal to the area where the decreased amplitude is noted. Also, measurements of arterial pressures are performed at each cuff site. A difference in pressure of greater than 20 mm Hg indicates a degree of arterial occlusion in the extremity. A positive result is reliable evidence of arteriosclerotic peripheral vascular occlusion. However, a negative result does not definitely exclude this diagnosis because extensive vascular collateralization can compensate for even a complete arterial occlusion.

An ankle/brachial ratio of less than 0.9 indicates peripheral vascular disease in the lower extremity. Arterial plethysmography can also be performed immediately after exercise to determine whether symptoms of claudication are caused by peripheral vascular occlusive disease.

Although it is not as accurate as arteriography (see p. 114), plethysmography is performed without serious complications and can be done for extremely ill patients who cannot be transported to the arteriography laboratory.

## 714 plethysmography, arterial

### Interfering factors

- Arterial occlusion proximal to the extremity
- Cigarette smoking can cause transient arterial constriction.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting is required.
- PT** Inform the patient that this test is painless.
- PT** Tell the patient that he or she must lie still during the testing procedure.
  - Remove all clothing from the patient's extremities.
- PT** Instruct the patient to avoid smoking for at least 30 minutes before the test. Nicotine creates constriction of the peripheral arteries and alters the test results.

#### During

- Note the following procedural steps:
  1. The patient is placed in the semirecumbent position.
  2. The cuffs are applied to the extremities and then inflated to 65 mm Hg to increase their sensitivity to pulse waves.
  3. The pulse waves are recorded on plethysmographic paper.
  4. The amplitudes and form of the pulse wave of each cuff are measured and compared. A marked reduction in wave amplitude indicates arterial occlusive disease.
- This noninvasive test is usually performed in the vascular laboratory or at the patient's bedside by a vascular technologist in approximately 30 minutes.
- PT** Inform the patient that results are usually interpreted by a physician and are available in a few hours.
- PT** Tell the patient that no discomfort is associated with testing.

#### After

- Encourage the patient to verbalize concerns regarding results.

### Abnormal findings

Arterial occlusive disease

Arterial trauma

Small vessel diabetic changes

Vascular diseases (e.g., Raynaud phenomenon)

Arterial embolization

#### notes

## pleural biopsy

**Type of test** Microscopic examination of tissue

### Normal findings

No evidence of pathology

### Test explanation and related physiology

This test may be indicated when the pleural fluid obtained by thoracentesis (see p. 864) is exudative fluid, which suggests infection, neoplasm, or tuberculosis. The pleural biopsy can distinguish among these disease processes. It is also performed when chest imaging indicates a pleural-based tumor, reaction, or thickening.

*Pleural biopsy* is the removal of pleural tissue for histologic examination. It is usually performed by a percutaneous needle biopsy. It also can be performed via thoracoscopy (see p. 870). Pleural tissue also may be obtained by an *open pleural biopsy*, which involves a limited thoracotomy and requires general anesthesia. For this procedure, a small intercostal incision is made, and the biopsy of the pleura is done under direct observation. The advantage of an open procedure is that a larger piece of pleura can be obtained.

### Contraindications

- Patients with prolonged bleeding or clotting times

### Potential complications

- Bleeding or injury to the lung
- Pneumothorax

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Obtain informed consent for this procedure.

**PT** Tell the patient that no fasting or sedation is required.

**PT** Instruct the patient to remain very still during the procedure. Any movement may cause inadvertent needle damage.

#### During

- Note the following procedural steps for percutaneous needle biopsy:
  1. This procedure is usually performed with the patient in a sitting position with his or her shoulders and arms elevated and supported by a padded overbed table.

## 716 pleural biopsy

2. After the presence of fluid has been determined by the thoracentesis technique, the skin overlying the biopsy site is anesthetized and pierced with a scalpel blade.
  3. A needle is inserted with a cannula until fluid is removed. (Some fluid is left in the pleural space after the thoracentesis to make the biopsy easier.)
  4. The inner needle is removed, and a blunt-tipped, hooked biopsy trocar attached to a three-way stopcock is inserted into the cannula.
  5. The patient is instructed to expire all air and then perform the Valsalva maneuver to prevent air from entering the pleural space.
  6. The cannula and biopsy trocar are withdrawn while the hook catches the parietal wall and takes a specimen with its cutting edge.
  7. Usually three biopsy specimens are taken from different sites during the same session.
  8. The specimens are placed in a fixative solution and sent to the laboratory immediately.
  9. After the specimens are taken, additional parietal fluid can be removed.
- Note that this procedure is performed by a physician at the patient's bedside, in a special procedure room, or in the physician's office in approximately 30 minutes.
- PT** Tell the patient that because of the local anesthetic little discomfort is associated with this procedure.

### After

- Apply an adhesive bandage to the biopsy site.
  - Note that a chest x-ray image is usually taken to detect the potential complication of pneumothorax.
  - Observe the patient for signs of respiratory distress (e.g., diminished breath sounds) on the side of the biopsy.
- PT** Instruct the patient to report any shortness of breath.
- Observe the patient's vital signs frequently for evidence of bleeding (increased pulse, decreased blood pressure).
  - Ensure that the biopsy specimen is sent immediately to the laboratory.

### Abnormal findings

Neoplasm  
Tuberculosis

### notes



## porphyrins and porphobilinogens

**Type of test** Urine (fresh and 24 hour)

### Normal findings

	Male (mcg/24 hr)	Female (mcg/24 hr)
Total porphyrins	8-149	3-78
Uroporphyrin	4-46	3-22
Coproporphyrin	< 96	< 60

Porphobilinogens: 0-2 mg/24 hr or 0-8.8  $\mu$ mol/day (SI units)

### Test explanation and related physiology

This test is a quantitative measurement of porphyrins and porphobilinogens. Along with measurement of aminolevulinic acid (see p. 321), the various forms of porphyria can be identified.

Porphyria is a group of genetic disorders associated with enzyme deficiencies involved with porphyrin synthesis or metabolism. Porphyrins (e.g., uroporphyrin and coproporphyrin) and porphobilinogens are important building blocks in the synthesis of heme. Heme is incorporated into hemoglobin in the erythroid cells. In most forms of porphyria, increased levels of porphyrins and porphobilinogens are found in the urine. Heavy metal (lead) intoxication is also associated with increased porphyrins in the urine.

Urine tests for porphyrins are not as accurate as plasma measurements and pattern identification for the various forms of porphyria. They are accurate, however, in screening for porphyria, especially the intermittent variety. Porphyrin fractionation of erythrocytes and plasma provides specific assays for primary RBC porphyrins. These assays are predominantly used to differentiate the various forms of congenital porphyrias. Plasma measurement of *free erythrocyte protoporphyrin (FEP)* is helpful in the diagnosis of iron-deficiency anemia or lead intoxication.

### Interfering factors

- ✚ Drugs that may alter test results include aminosalicic acid, barbiturates, chloral hydrate, chlorpropamide, ethyl alcohol, griseofulvin, morphine, oral contraceptives, phenazopyridine, procaine, and sulfonamides.

## 718 porphyrins and porphobilinogens

### Procedure and patient care

- See inside front cover for Routine Urine Testing for random and 24-hour collection.
- Protect the specimen from light.
- PT** Instruct the patient to avoid alcohol during the collection period.
- Keep the 24-hour urine in a light-resistant specimen bottle with a preservative to prevent degradation of the light-sensitive porphyrin.
- PT** Encourage the patient to drink fluids during the 24 hours unless contraindicated for medical purposes.

### Abnormal findings

#### ▲ Increased levels

Porphyrias  
Liver disease  
Lead poisoning  
Pellagra

#### notes

## positron emission tomography (PET scan)

**Type of test** Nuclear scan, x-ray

### Normal findings

No abnormal areas of increased or decreased uptake

### Test explanation and related physiology

PET scanning is used in many areas of medicine. The greatest use of PET scan has been in the fields of neurology, cardiology, and oncology.

In PET scanning, radioactive chemicals are administered to the patient. These chemicals are used in the normal metabolic process of the cells of the particular organ being imaged. Positrons emitted from the radioactive chemicals in the organ are sensed by a series of detectors positioned around the patient. The positron emissions are recorded and reconstructed by computer analysis into a high-resolution three-dimensional image indicating a particular metabolic process in a specific anatomic site. A CT x-ray scan (see p. 278) is performed on the patient at the same time to assist in the development and interpretation of the images created. PET/CT scans provide images representing not only anatomy but also physiology.

Depending on the particular radionuclide used, PET can demonstrate the glucose metabolism, oxygenation, blood flow, and tissue perfusion of any specific area. Pathologic conditions are recognized and diagnosed by alterations in the normal metabolic process.

### Neurology

Pathologic areas of the brain that are more metabolically active (e.g., cancers) more avidly take up fluorodeoxyglucose (FDG) than do normal areas. Because of the high physiologic rate at which glucose is metabolized by normal brain tissue, the detectability of tumors with only modest increases in glucose metabolism, such as low-grade tumors and, in some cases, recurrent tumors, is difficult with FDG.

Epilepsy, PD, and Huntington disease are identified as localized areas of increased metabolic activity indicating rapid nerve firing. Brain trauma resulting in a hematoma or bleeding is evident as decreased metabolic activity in the area of trauma. Stroke can also be identified and its extent determined. With the use of radioactive water, brain blood flow can be determined. Areas of

## 720 positron emission tomography

decreased blood flow take up less  $H_2O$  than normal areas and represent areas at risk for stroke.

Alzheimer disease can be recognized by identifying hypometabolism in multiple areas of the brain (temporal and parietal lobe). PET scanning with amyloid imaging has been very helpful in identifying amyloid protein precursors (see p. 58) in the brain in patients with Alzheimer disease. A negative PET scan with amyloid imaging eliminates the possibility of Alzheimer disease in a patient with cognitive impairment. Because other neurologic conditions (especially in elderly people) are also associated with amyloid neuritic plaques, a positive scan does not establish the diagnosis of Alzheimer disease with certainty.

### Cardiology

PET scans of the heart are one of the best tests to indicate decreased blood flow, indicating coronary artery occlusive disease. PET scans are also used when cardiac muscle function is reduced. A PET scan can indicate whether the dysfunction arises from reversible ischemic muscle that would benefit from revascularization or from muscle tissue that is no longer viable. In the former case, surgical revascularization should be considered. In the latter case, revascularization would not be beneficial.

### Oncology

The most commonly used agent in oncology is FDG because increased glucose metabolism is so prevalent in malignant tumors compared with normal or benign pathologic tissue. PET can be used to visualize rapidly growing tumors and indicate their anatomic location. It is used to determine tumor response to therapy, identify recurrence of tumor after surgical removal, and differentiate tumor from other pathologic conditions (e.g., infection). PET is particularly helpful in identifying regional and metastatic spread for a particular tumor. In lung cancer, for example, if the FDG fails to concentrate in any area other than the primary tumor, no spread is suspected and the patient is considered an ideal candidate for surgery. PET has also been particularly useful for identifying metastasis from lung, melanoma, breast, pancreas, colon, lymphoma, and brain cancers.

Rapidly growing tumors are associated with a high metabolic rate and therefore concentrate FDG particularly well. The amount of uptake of FDG is measured by the standardized uptake value (SUV)—the amount of uptake of FDG in tumor compared with the normal tissue in that same area. SUV helps to

distinguish between benign and malignant lesions—the higher the SUV, the more likely the tumor is malignant.

When the SUV is greater than the “cutoff value” (as determined by each institution), cancer rather than a benign pathologic condition is suspected. PET scanning is particularly helpful in the evaluation of solitary pulmonary nodules.

## Bone

A PET/CT scan with a sodium fluoride F18 injection ( $^{18}\text{F}$  NaF) scans the entire skeletal system and produces high-resolution images of the bones. These images are used to detect areas of abnormal bone growth associated with tumors. This test is more accurate than conventional nuclear bone scans. The PET/CT scan of the bone is particularly helpful for patients with prostate or breast cancer. The uptake of  $^{18}\text{F}$  NaF in the skeleton reflects sites of increased blood flow and bone remodeling associated with bone injury or metastatic disease.

Small parts PET scans are being used for foot inflammatory pathology. *PET mammography or positron emission mammography (PEM)* is seeing growing use as a tool for diagnostic breast imaging. PEM holds the promise of improving the sensitivity and specificity of routine mammography (see p. 606).

## Interfering factors

- Recent use (within 24 hours) of caffeine, alcohol, or tobacco may affect test results.
- Ingestion of a small- to moderate-sized meal can cause a marked uptake of FDG in the gut and muscles, thereby leaving little or no radionuclide to be taken up by tumor. This can cause a false-negative result.
- Anxiety can cause increased uptake in multiple areas (e.g., neck, upper mediastinum) of the body. If the patient is anxious, sedatives can be administered 30 minutes before testing. However, these could interfere with PET scanning of the brain if cognitive activities will be used to measure changes in brain activity.
- Mild to moderate exercise can instigate marked uptake of FDG in the muscles, thereby leaving little or no radionuclide to be taken up by tumor. This causes a false-negative result.
- The liver and spleen avidly take up FDG. Therefore these organs are difficult to evaluate on PET imaging.
- FDG is excreted by the urinary system. As a result, the bladder may obscure areas of increased uptake in the pelvis.

## 722 positron emission tomography

- Uptake of FDG can occur in the lymph node basin draining the site of the FDG injection. If PET is being done to stage tumors that could metastasize to those lymph nodes, the FDG should be injected on the contralateral side.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- Obtain informed consent if required by the institution.
- PT** Inform the patient that he or she may have an IV line inserted.
- PT** Inform the patient that he or she may need to restrict food or fluids for 4 hours on the day of the test. The patient should refrain from alcohol, caffeine, and tobacco for 24 hours.
- PT** Instruct patients with diabetes to take their pretest dose of insulin at a meal 3 to 4 hours before the test.
- PT** Tell the patient that no sedatives or tranquilizers should be taken, because he or she may need to perform certain mental activities during the brain PET scan.
- PT** Tell the patient to empty the bladder before the test for comfort. A Foley catheter may be inserted for PET scanning of the pelvic region.
- PT** Tell the patient that the only discomfort associated with this study is insertion of the IV line.
- Depending on the organ being evaluated, specific protocols may exist for the examination.

#### During

- Note the following procedural steps:
  1. The patient is positioned in a comfortable, reclining chair.
  2. The radioactive material can be infused through an IV line.
  3. The gamma rays that penetrate the tissues are recorded outside the body by a circular array of detectors and are displayed by a computer.
  4. Extraneous auditory and visual stimuli are minimized by a blindfold and ear plugs.
  5. If the chest is being scanned, instruct the patient to breathe in a shallow manner until the middle of the chest is reached. Then ask the patient to hold the breath after expiration until the middle of the abdomen is reached. This will improve visibility of the chest anatomy.
- Note that a physician performs this procedure with a trained technologist in approximately 40 to 90 minutes.

**After**

- PT** Instruct the patient to change position slowly from lying to standing to avoid postural hypotension.
- PT** Encourage the patient to drink fluids and urinate frequently to aid in removal of the radioisotope from the bladder.

**Abnormal findings**

Myocardial infarction	Huntington disease
Coronary artery disease	Dementia
Cerebrovascular accident	Alzheimer disease
Epilepsy	Malignant tumor
PD	

**notes**

**potassium, blood (K)****Type of test** Blood**Normal findings**

Adult/elderly: 3.5-5.0 mEq/L or 3.5-5.0 mmol/L (SI units)

Child: 3.4-4.7 mEq/L

Infant: 4.1-5.3 mEq/L

Newborn: 3.9-5.9 mEq/L

**Possible critical values**

Adult: &lt; 2.5 or &gt; 6.5 mEq/L

Newborn: &lt; 2.5 or &gt; 8 mEq/L

**Test explanation and related physiology**

Potassium (K) is the major cation within the cell. Normal serum K concentration is approximately 4 mEq/L. Because the serum concentration of K is so small, minor changes in concentration have significant consequences. K is excreted by the kidneys, and there is no resorption of K from the kidneys. Therefore if K is not adequately supplied in the diet (or by IV administration in patients who are unable to eat), serum K levels can drop rapidly.

Serum K concentration depends on many factors, including the following:

- *Aldosterone*. This hormone tends to increase renal losses of K.
- *Sodium resorption*. As sodium is resorbed, K is lost.
- *Acid-base balance*. Alkalotic states tend to lower serum K levels by causing a shift of K into the cell. Acidotic states tend to raise serum K levels by reversing that shift.

An electrocardiogram may demonstrate peaked T waves, a widened QRS complex, and depressed ST segment in hyperkalemia. Hypokalemia is associated with increased cardiac sensitivity to digoxin, cardiac arrhythmias, flattened T waves, and prominent U waves. This electrolyte has profound effects on the heart rate and contractility. The K level should be followed carefully in patients with uremia, Addison disease, vomiting, or diarrhea; in patients on steroid therapy; and in patients taking K-depleting diuretics. K must be closely monitored in patients taking digitalis-like drugs because cardiac arrhythmias may be induced by hypokalemia and digoxin.



## Interfering factors

- Movement of the forearm with a tourniquet in place may increase K levels.
- Hemolysis during venipuncture causes increased levels.
- ✚ Drugs that may cause *increased* K levels include aminocaproic acid, antibiotics, antineoplastic drugs, captopril, epinephrine, heparin, histamine, isoniazid, lithium, mannitol, potassium-sparing diuretics, potassium supplements, and succinylcholine.
- ✚ Drugs that may cause *decreased* K levels include acetazolamide, aminosalicic acid, amphotericin B, carbenicillin, cisplatin, diuretics, glucose infusions, insulin, laxatives, lithium carbonate, penicillin G sodium (high doses), phenothiazines, salicylates, and sodium polystyrene sulfonate.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: no
  - Blood tube commonly used: red or green
- PT** Instruct the patient to avoid opening and closing the hand after a tourniquet is applied.
- Avoid hemolysis.
  - Evaluate patients with increased or decreased K levels for cardiac arrhythmias.
  - Monitor patients taking digoxin and diuretics for hypokalemia.
  - If indicated, administer resin exchanges (e.g., Kayexalate enema) to correct hyperkalemia.

## Abnormal findings

### ▲ Increased levels (hyperkalemia)

- Excessive dietary intake
- Excessive IV intake
- Acute or chronic renal failure
- Hypoaldosteronism
- Aldosterone-inhibiting diuretics
- Crush injury to tissues
- Hemolysis
- Transfusion of hemolyzed blood
- Infection
- Acidosis
- Dehydration

### ▼ Decreased levels (hypokalemia)

- Deficient dietary intake
- Deficient IV intake
- Burns
- GI disorders (e.g., diarrhea)
- Diuretics
- Hyperaldosteronism
- Cushing syndrome
- Renal tubular acidosis
- Licorice ingestion
- Insulin administration
- Glucose administration
- Ascites
- Renal artery stenosis
- Cystic fibrosis
- Trauma
- Surgery

notes

## potassium, urine (K)

**Type of test** Urine (24-hour)

### Normal findings

25-100 mEq/L/day or 25-100 mmol/day (SI units)  
(Values vary greatly with diet.)

### Test explanation and related physiology

K is the major cation within the cell (see p. 724). K can be measured in both a spot collection and a 24-hour urine collection. K concentration in the urine depends on many factors. Aldosterone and, to a lesser extent, glucocorticosteroids tend to increase the renal losses of K. Acid-base balance is dependent on K excretion to a small degree. In alkalotic states, hydrogen can be resorbed in exchange for K. The kidneys cannot resorb K. Therefore the intake of K is balanced by kidney excretion.

### Interfering factors

- Dietary intake affects K levels.
- Excessive intake of licorice may cause increased levels.
- ✖ Drugs that may cause *increased* levels include diuretics, glucocorticoids, and salicylates.

### Procedure and patient care

- See inside front cover for Routine Urine Testing.

### Abnormal findings

#### ▲ Increased levels

Chronic renal failure  
Renal tubular acidosis  
Cushing syndrome  
Hyperaldosteronism  
Excessive intake of licorice  
Alkalosis  
Diuretic therapy

#### ▼ Decreased levels

Dehydration  
Addison disease  
Malnutrition  
Vomiting  
Diarrhea  
Malabsorption  
Acute renal failure

notes

**prealbumin** (PAB, Thyroxine-binding prealbumin [TBPA], Thyretin, Transthyretin)

**Type of test** Blood; urine (24-hour); cerebrospinal fluid (CSF) analysis

### Normal findings

#### Serum

Adult/elderly: 15-36 mg/dL or 150-360 mg/L (SI units)

#### Child

< 5 days: 6-21 mg/dL

1-5 years: 14-30 mg/dL

6-9 years: 15-33 mg/dL

10-13 years: 22-36 mg/dL

14-19 years: 22-45 mg/dL

#### Urine (24-hour)

0.017-0.047 mg/day

#### CSF

Approximately 2% of total CSF protein

### Possible critical values

Serum prealbumin levels less than 10.7 mg/dL indicate severe nutritional deficiency.

### Test explanation and related physiology

Prealbumin is one of the major plasma proteins. Because prealbumin can bind thyroxine, it is also called *thyroxine-binding prealbumin* (TBA). However, prealbumin is secondary to thyroxine-binding globulin in the transportation of triiodothyronine ( $T_3$ ) and thyroxine ( $T_4$ ). Prealbumin also plays a role in the transport and metabolism of vitamin A.

Because prealbumin levels in serum fluctuate more rapidly in response to alterations in synthetic rate than do those of other serum proteins, clinical interest in the quantification of serum prealbumin has centered on its usefulness as a marker of nutritional status. Its half-life of 1.9 days is much less than the 21-day half-life of albumin (see p. 746). Because prealbumin has a short half-life, it is a sensitive indicator of any change affecting protein synthesis and catabolism. For this reason, prealbumin is frequently ordered to monitor the effectiveness of total parenteral nutrition (TPN).

Prealbumin is significantly reduced in hepatobiliary disease because of impaired synthesis. Prealbumin is also a negative acute-phase reactant protein; serum levels decrease in inflammation,

malignancy, and protein-wasting diseases of the intestines or kidneys. Because zinc is required for synthesis of prealbumin, low levels occur with a zinc deficiency. Increased levels of prealbumin occur in Hodgkin disease and chronic kidney disease.

### Interfering factors

- Coexistent inflammation may make the interpretation of test results impossible.
- ▼ Drugs that may cause *increased* levels include anabolic steroids, androgens, and prednisolone.
- ▼ Drugs that may cause *decreased* levels include amiodarone, estrogens, and oral contraceptives.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- If the patient is going to collect a 24-hour urine specimen, provide a collection bottle. (See guidelines on inside front cover.)

### Abnormal findings

#### ▲ Increased levels

Some cases of nephrotic syndrome  
Hodgkin disease  
Chronic kidney disease  
Pregnancy

#### ▼ Decreased levels

Malnutrition  
Liver damage  
Burns  
Salicylate poisoning  
Inflammation  
Infection

notes

## pregnancy-associated plasma protein-A (PAPP-A)

**Type of test** Blood

### Normal findings

Normal values vary by laboratory and duration of pregnancy

### Test explanation and related physiology

*Pregnancy-associated plasma protein-A (PAPP-A)* is made by the trophoblasts and released into the maternal circulation during pregnancy. Women with low blood levels of PAPP-A at 8 to 14 weeks of gestation have an increased risk of intrauterine growth restriction, trisomy 21, premature delivery, preeclampsia, and stillbirth. This protein rapidly rises in the first trimester of normal pregnancy. However, in Down syndrome-affected pregnancy, serum levels are half that of unaffected pregnancies. Furthermore, low levels of PAPP-A in maternal serum in the first trimester are associated with adverse fetal outcomes, including fetal death in utero and intrauterine growth retardation. This test is commonly used with other pregnancy and screening tests (see p. 612).

PAPP-A is present in unstable atherosclerotic plaques, and circulating levels are elevated in acute coronary syndromes, which may reflect the instability of the plaques. PAPP-A is an independent marker of unstable angina and acute myocardial infarction.

### Interfering factors

- Levels increase with increased maternal body weight and longer duration of pregnancy.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Allow the patient to express her concerns and fears regarding the potential for birth defects.

**PT** Assist the patient in scheduling and obtaining more accurate diagnostic testing if the results are positive.

### Abnormal findings

Positive screening tests (trisomy 21, trisomy 18, neural tube defects, abdominal wall defects)

Coronary atherosclerotic disease

### notes

## pregnanediol

**Type of test** Urine (24-hour)

### Normal findings

< 2 years: < 0.1 mg/day

< 9 years: < 0.5 mg/day

10-15 years: 0.1-1.2 mg/day

Adult male: 0-1.9 mg/day

Adult female

Follicular phase: < 2.6 mg/day

Luteal: 2.6-10.6 mg/day

Pregnancy

First trimester: 10-35 mg/day

Second trimester: 35-70 mg/day

Third trimester: 70-100 mg/day

### Test explanation and related physiology

Urinary pregnanediol is measured to evaluate progesterone production by the ovaries and placenta. Both serum progesterone levels and the urine concentration of progesterone metabolites (pregnanediol and others) are significantly increased during the latter half of an ovulatory cycle. Pregnanediol is the most easily measured metabolite of progesterone.

Because pregnanediol levels rise rapidly after ovulation, this study is useful in documenting whether ovulation has occurred and, if so, its exact time. During pregnancy, pregnanediol levels normally rise because of the placental production of progesterone. Repeated assays can be used to monitor the status of the placenta in women who are having difficulty becoming pregnant or maintaining a pregnancy. This study is also used to monitor high-risk pregnancies.

Hormone assays for urinary pregnanediol are primarily used today to monitor progesterone supplementation in patients with an inadequate luteal phase. Urinary assays may be supplemented by plasma assays (progesterone assay; see p. 733), which are quicker and more accurate.

### Interfering factors

- ✦ Drugs that may cause *increased* levels include adrenocorticotrophic hormone.
- ✦ Drugs that may cause *decreased* levels include oral contraceptives and progesterones.

### **Procedure and patient care**

- See inside front cover for Routine Urine Testing.
- Record on the laboratory slip the date of the last menstrual period or the week of gestation during pregnancy.

### **Abnormal findings**

#### **▲ Increased levels**

Ovulation  
Pregnancy  
Luteal cysts of ovary  
Arrhenoblastoma of ovary  
Hyperadrenocorticalism  
Choriocarcinoma of ovary  
Adrenocortical hyperplasia

#### **▼ Decreased levels**

Threatened abortion  
Fetal death  
Toxemia of pregnancy  
Amenorrhea  
Ovarian hypofunction  
Placental failure  
Preeclampsia  
Ovarian neoplasm  
Breast neoplasm

**notes**



## progesterone assay

**Type of test** Blood

### Normal findings

< 9 years: < 20 ng/dL

10-15 years: < 20 ng/dL

Adult male: 10-50 ng/dL

Adult female

Follicular phase: < 50 ng/dL

Luteal: 300-2500 ng/dL

Postmenopausal: < 40 ng/dL

Pregnancy

First trimester: 725-4400 ng/dL

Second trimester: 1950-8250 ng/dL

Third trimester: 6500-22,900 ng/dL

### Test explanation and related physiology

The major effect of progesterone is to induce the development of the secretory phase of the endometrium in anticipation of implantation of a fertilized ovum. Normally, progesterone is secreted by the ovarian corpus luteum after ovulation. Serum progesterone level is significantly increased during the second half of the ovulatory cycle. Normally, blood samples drawn at days 8 and 21 of the menstrual cycle show a large increase in progesterone levels in the latter specimen, indicating that ovulation has occurred. Therefore this study is useful in documenting whether ovulation has occurred and, if so, its exact time. This is very useful information for a woman who has difficulty becoming pregnant.

In pregnancy, progesterone is produced by the corpus luteum for the first few weeks. After that, the placenta begins to make progesterone. Progesterone levels should progressively rise during pregnancy because of placental production. Repeated assays can be used to monitor the status of the placenta in cases of high-risk pregnancy. Progesterone assay is also used today to monitor progesterone supplementation in patients with an inadequate luteal phase to maintain an early pregnancy.

### Interfering factors

- Recent use of radioisotopes may affect test results.
- Hemolysis caused by rough handling of the sample may affect test results.

## 734 progesterone assay

- ▼ Drugs that may interfere with test results include clomiphene, estrogen, and progesterone.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Indicate the date of the last menstrual period on the laboratory slip.

### Abnormal findings

#### ▲ Increased levels

Ovulation  
Pregnancy  
Luteal cysts of ovary  
Hyperadrenocorticalism  
Adrenocortical hyperplasia  
Choriocarcinoma of ovary  
Hydatidiform mole of the uterus

#### ▼ Decreased levels

Preeclampsia  
Toxemia of pregnancy  
Threatened abortion  
Placental failure  
Fetal death  
Ovarian neoplasm  
Amenorrhea  
Ovarian hypofunction

### notes

## progesterone receptor assay (PR assay, PRA, PgR)

**Type of test** Tumor-specimen analysis

### Normal findings

Immunochemistry

Negative: < 5% of the cells stain for receptors

Positive: > 5% of the cells stain for receptors

Reverse-transcriptase polymerase chain reaction (RT-PCR)

Negative: < 5.5 units

Positive: > 5.5 units

### Test explanation and related physiology

The PR assay is used in determination of the prognosis and treatment of breast cancer and, to a lesser degree, other cancers. This assay helps determine whether a tumor is likely to respond to endocrine medical or surgical therapy. The test is done on breast cancer specimens when a primary or recurrent cancer is identified; it is usually done in conjunction with an estrogen receptor (ER) assay (see p. 389) to increase the predictability of a tumor response to hormone therapy. Tumor response rates to medical or surgical hormonal manipulation are potentiated if the ER assay is positive. Tumor response rates are as follows:

- ER positive, PR positive: 75%
- ER negative, PR positive: 60%
- ER positive, PR negative: 35%
- ER negative, PR negative: 25%

Other tumors (e.g., ovarian, melanoma, uterine, pancreatic) are occasionally studied for ER and PR assay. This is done mostly within clinical trials.

### Interfering factors

- ✦ The use of such hormones as progesterone or estrogen may cause false-negative results.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Prepare the patient for breast biopsy per routine protocol.
- Record the menstrual status of the patient.
- Record any exogenous hormone the patient may have used during the past 2 months.

## 736 progesterone receptor assay

**PT** Instruct the patient to discontinue exogenous hormones before breast biopsy. This is done in consultation with the physician.

### During

- The surgeon obtains tissue.
- This tissue should be placed on ice or in formalin.
- Part of the tissue is used for routine histology. A portion of the paraffin block is sent to a reference laboratory.
- Results are usually available in 1 week.

### After

- Provide routine postoperative care.

### Abnormal findings

#### ▲ Increased levels

Hormonally dependent cancer

notes

## prolactin levels (PRLs)

**Type of test** Blood

### Normal findings

Adult male: 3-13 ng/mL

Adult female: 3-27 ng/mL

Pregnant female: 20-400 ng/mL

### Test explanation and related physiology

Prolactin is a hormone secreted by the anterior pituitary gland (adenohypophysis). In females, prolactin promotes lactation. Its role in males is not clear. During sleep, prolactin levels increase two- to threefold to circulating levels equaling those of pregnant women. With breast stimulation, pregnancy, nursing, stress, or exercise, a surge of this hormone occurs. It is elevated in patients with prolactin-secreting pituitary acidophilic or chromophobic adenomas. To a lesser extent, moderately high prolactin levels have been observed in women with secondary amenorrhea (i.e., postpubertal) and galactorrhea. Paraneoplastic tumors (e.g., lung cancer) may cause ectopic secretion of prolactin as well. In general, very high prolactin levels are more likely to be caused by pituitary adenoma than other causes.

The prolactin level is helpful for monitoring the disease activity of pituitary adenomas. Several *prolactin stimulation tests* (with thyrotropin-releasing hormone or chlorpromazine) and *prolactin suppression tests* (with levodopa) have been designed to help differentiate pituitary adenoma from other causes of prolactin overproduction.

### Interfering factors

- Stress from illness, trauma, surgery, or even the fear of a blood test can elevate levels.
- Prolactin levels are transiently elevated after seizures.
- ✚ Drugs that may cause *increased* values include anticonvulsants, antihistamines, anti-nausea/antiemetic drugs, antipsychotic drugs, antituberculosis medications, ergot derivatives, estrogens/progesterone, histamine antagonists, monoamine oxidase inhibitors, opiates, oral contraceptives, reserpine, serotonin reuptake inhibitors, several antihypertensive drugs, and some illegal drugs.
- ✚ Drugs that may cause *decreased* values are clonidine, dopamine, ergot alkaloid derivatives, and levodopa.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- ▶ Inform the patient that blood should be drawn in the morning.
- Transfer the specimen to the laboratory as soon as possible. If a delay occurs, the specimen should be placed on ice.

## Abnormal findings

### ▲ Increased levels

Galactorrhea  
Amenorrhea  
Prolactin-secreting pituitary tumor  
Infiltrative diseases of the hypothalamus and pituitary stalk  
Hypothyroidism  
Renal failure  
Anorexia nervosa  
Perineoplastic ectopic production of prolactin  
Metastatic cancer to the pituitary gland  
Polycystic ovary syndrome  
Stress  
Empty sella syndrome

### ▼ Decreased levels

Pituitary apoplexy  
Pituitary destruction from tumor (craniopharyngioma)

## notes

## ProstaScint scan (Radioimmunosциntigraphy [RIS])

**Type of test** Nuclear scan

### Normal findings

No increased uptake outside the prostate gland

### Test explanation and related physiology

By using a radionuclide that is able to attach to prostate cancer cells only, metastatic prostate cancer outside the prostate gland can be easily identified. Disease outside the prostate (e.g., retroperitoneum, liver, lung, or bone) indicates metastatic prostate cancer.

This scan is helpful in staging newly diagnosed prostate cancer patients who are at high risk for metastatic disease to the lymph nodes or other organs. This test can also be used to identify recurrent or metastatic disease after curative therapy. This test can document completeness of anti-prostate cancer therapy.

### Contraindications

- Patients who are unable to be immobile for 1 hour

### Interfering factors

- Prior bone scans confound image interpretation.
- Areas of inflammation (e.g., degenerative joint disease, inflammatory bowel disease, or recent trauma) can confound interpretation.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Explain that no fasting is required before the test.
- PT** Instruct the patient to use a mild laxative the evening before imaging. A cleansing enema may be administered 1 hour before imaging. Early in the scanning, radionuclide accumulates in the bowel.
- PT** Tell the patient to void before each image.

#### During

- After proper identification, the patient is injected with the radiolabeled monoclonal antibody.
- Initial images are obtained 30 minutes after injection. Images are repeated over as many as 5 days.

## 740 ProstaScint scan

- Little or no discomfort is associated with this procedure.
- The procedure takes approximately 1 hour each day over a period of 1 to 5 days.
- This procedure is performed in the nuclear medicine department.

### After

- PT Inform the patient that no precautions need to be taken by others against radiation exposure because only tracer doses of radioisotopes are used.
- Encourage the patient to increase oral fluid intake on testing days.

### Abnormal findings

Primary or recurrent prostate cancer

notes



## prostate/rectal sonogram (Ultrasound prostate)


**Type of test** Ultrasound

### Normal findings

Normal size, contour, and consistency of the prostate gland

### Test explanation and related physiology

Rectal ultrasonography of the prostate is a very valuable tool in the early diagnosis of prostate cancer. When combined with rectal digital examination and prostate-specific antigen (see p. 743), very small prostate cancers can be identified. Prostate/rectal sonography is also helpful in evaluating the seminal vessels and other perirectal tissue. Ultrasound is very helpful in guiding the direction of a prostate biopsy (Figure 34) and in quantitating the volume of prostate cancer. When radiation therapy implantation is required for treatment, ultrasonography is used to map the exact location of the prostate cancer. Rectal ultrasonography is very helpful in staging rectal cancers as well. The depth of

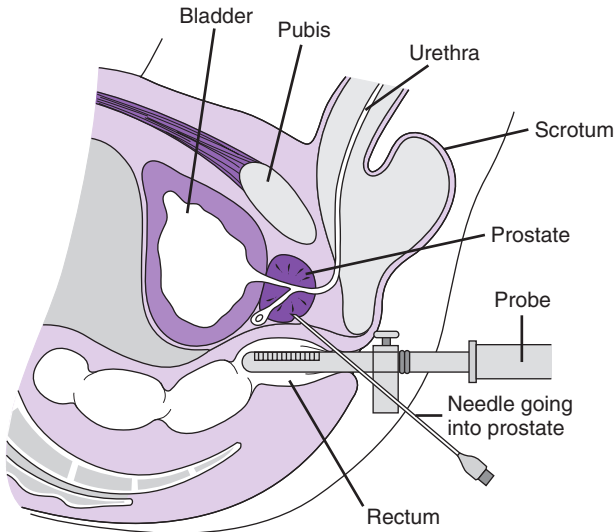

**FIG. 34** Rectal ultrasonography. Diagram demonstrating transrectal biopsy of the prostate.

## 742 prostate/rectal sonogram

transmural involvement and presence of extrarectal extension can be accurately assessed.

### Contraindications

- Patients with latex allergy. Rectal ultrasonography requires placement of the probe in a latex condom–like sac.

### Interfering factors

- Stool within the rectum

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

**PT** Instruct the patient that a small-volume rectal enema may be required approximately 1 hour before the ultrasound examination.

#### During

- The patient is placed in the left lateral decubitus position.
- A digital rectal examination may be performed to assess the prostate gland or rectal tumor.
- A draped and lubricated ultrasound probe is placed in the rectum.
- Scans are performed in various spatial planes.

#### After

- Provide the patient with tissue material to cleanse the perianal area.

### Abnormal findings

Prostate cancer

Benign prostatic hypertrophy

Prostatitis

Seminal vesicle tumor

Prostate abscess

Perirectal abscess

Intrarectal or perirectal tumor

### notes

**prostate-specific antigen (PSA), prostate cancer gene 3 (PCA3), MI-prostate score (MiPS)**

**Type of test** Blood

**Normal findings**

0-2.5 ng/mL is low

2.6-10 ng/mL is slightly to moderately elevated

10-19.9 ng/mL is moderately elevated

≥ 20 ng/mL is significantly elevated

**Test explanation and related physiology**

PSA is elevated in patients with prostate diseases (e.g., cancer, infection, and benign hypertrophy). PSA can be detected in all men; however, levels can be greatly increased in patients with prostate cancer. The higher the levels, the greater the tumor burden. The PSA assay is also a sensitive test for monitoring response to therapy. Significant elevation in PSA subsequently indicates the recurrence of prostatic cancer.

There is considerable controversy regarding the use of PSA screening among asymptomatic men. Most professionals do not recommend the routine use of PSA screening in men of any age. Approximately 80% of PSA screening testing is falsely positive. A positive screening test result may trigger more aggressive and costly testing. It is important to be aware that some patients with early prostate cancer will not have elevated levels of PSA. PSA levels also may be minimally elevated in patients with benign prostatic hypertrophy (BPH) and prostatitis. In an effort to increase the accuracy of PSA testing, other measures of PSA or prostate proteins can be used:

- *PSA velocity*: PSA velocity is the change in PSA levels over time. A sharp rise in the PSA level raises the suspicion of cancer and may indicate a fast-growing cancer.
- *Age-adjusted PSA*: Age is an important factor in increasing PSA levels.
- *PSA density*: PSA density considers the relationship of the PSA level to the size of the prostate.
- *Free*: Free (protein unbound) PSA is a more accurate test for screening.
- *Prostate-specific proteins*: Prostate proteins (*prostatic-specific membrane antigen or early prostate cancer antigen [EPCA]*) may represent improved markers for prostate cancer.

## 744 prostate-specific antigen

- *Prostate cancer-specific biomarkers*: The most commonly tested marker is prostate cancer gene 3 (PCA3). Other genetic markers tested include GOLPH2, SPINK1, and TMPRSS2-ERG. These biomarkers are more specific to cancer and can indicate the risk that a prostate biopsy will show cancer.
- *PSA isoforms*: These isoforms can further increase the specificity of PSA to identify prostate cancer. The p2PSA is particularly specific.
- *Prostate Health Index (phi)*: This is particularly helpful for patients whose total PSA is between 4 and 10 ng/mL. This is a formula that is calculated as follows:

$$\text{phi} = \text{p2PSA} / \text{freePSA} \times \text{total PSA}^2.$$

Calculating phi can reduce the number of unnecessary prostate biopsies by 30%.

- *Mi-Prostate Score (MiPS)* combines the serum PSA with TMPRSS2-ERG and PCA3 in the urine. MiPS is used to predict a patient's risk for having prostate cancer detected by standard biopsy. The test also predicts the patient's risk for having potentially aggressive prostate cancer (Gleason score > 6).
- *Genomic prostate assay* (see p. 453, genomic testing), when used with commonly used pathology markers, can be used to indicate the likely clinical course of a known prostate cancer. *Oncotypedx* and *Promark* testing are commonly used tests. A low score is indicative of a less aggressive cancer for which observation alone may be reasonable. A high score indicates a more aggressive cancer that should be treated with aggressive intervention such as surgery or radiation.

### Interfering factors

- Rectal examinations may elevate PSA levels. The PSA specimen should be drawn before rectal examination of the prostate or several hours afterward.
- Prostatic manipulation by biopsy or transurethral resection of the prostate (TURP) may elevate PSA levels. The test should be done before surgery or 6 weeks afterward.
- Ejaculation within 24 hours of blood testing is associated with elevated PSA levels.
- Recent urinary tract infection or prostatitis can cause elevations of PSA for as long as 6 weeks.
- ✚ Diethylstilbestrol (DES) and finasteride (Propecia, Proscar) cause *decreased* levels of PSA by about 50%.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- The use of the percent-free PSA demands strict sample handling that is not required with the total PSA. Check for specific guidelines.

## Abnormal findings

### ▲ Increased levels

Prostate cancer

Benign prostatic hypertrophy

Prostatitis

notes

**protein** (Protein electrophoresis, Immunofixation electrophoresis [IFE], Serum protein electrophoresis [SPEP], Albumin, Globulin, Total protein)

**Type of test** Blood; urine

### Normal findings

Adult/elderly

Total protein: 6.4-8.3 g/dL or 64-83 g/L (SI units)

Albumin: 3.5-5 g/dL or 35-50 g/L (SI units)

Globulin: 2.3-3.4 g/dL

Alpha<sub>1</sub> globulin: 0.1-0.3 g/dL or 1-3 g/L (SI units)

Alpha<sub>2</sub> globulin: 0.6-1 g/dL or 6-10 g/L (SI units)

Beta globulin: 0.7-1.1 g/dL or 7-11 g/L (SI units)

Children

Total protein

Premature infant: 4.2-7.6 g/dL

Newborn: 4.6-7.4 g/dL

Infant: 6-6.7 g/dL

Child: 6.2-8 g/dL

Albumin

Premature infant: 3-4.2 g/dL

Newborn: 3.5-5.4 g/dL

Infant: 4.4-5.4 g/dL

Child: 4-5.9 g/dL

No protein abnormality on electrophoresis

### Test explanation and related physiology

Proteins are the most significant component contributing to the osmotic pressure in the vascular space. Albumin and globulin constitute most of the protein in the body and are measured together as the total protein.

*Albumin* makes up approximately 60% of the total protein. The major effect of albumin in the blood is to maintain colloidal osmotic pressure. Furthermore, albumin transports important blood constituents, such as drugs, hormones, and enzymes. Albumin is synthesized in the liver and is therefore a measure of hepatic function. When disease affects the liver cell, the hepatocyte loses its ability to synthesize albumin. The serum albumin level is greatly decreased. However, because the half-life of albumin is 12 to 18 days, severe impairment of

hepatic albumin synthesis may not be recognized until after that period.

*Globulins* represent all nonalbumin proteins. Alpha<sub>1</sub> globulins are mostly alpha<sub>1</sub> antitrypsin. Some transporting proteins, such as thyroid and cortisol-binding globulin, also contribute to this electrophoretic zone. Alpha<sub>2</sub> globulins include serum haptoglobins (which bind hemoglobin during hemolysis), ceruloplasmin (which is a carrier for copper), prothrombin, and cholinesterase (which is an enzyme used in the catabolism of acetylcholine). Beta<sub>1</sub> globulins include lipoproteins, transferrin, plasminogen, and complement proteins; beta<sub>2</sub> globulins include fibrinogen. Gamma globulins are the immunoglobulins (antibodies) (see p. 536). To a lesser degree, globulins also act as transport vehicles.

Serum albumin and some globulins are measures of nutrition. Malnourished patients, especially after surgery, have a greatly decreased level of serum proteins. Burn patients and patients who have protein-losing enteropathies and uropathies have low levels of protein despite normal synthesis.

In some diseases, albumin is selectively diminished, and globulins are normal or increased to maintain a normal total protein level. For example, in collagen vascular diseases (e.g., lupus erythematosus), capillary permeability is increased. Albumin, a molecule that is generally smaller than globulin, is selectively lost into the extravascular space. Another group of diseases similarly associated with low albumin, high globulin, and normal total protein levels is chronic liver diseases. In these diseases, the liver cannot produce albumin, but globulin is adequately made in the reticuloendothelial system. In both of these types of diseases, the albumin level is low, but the total protein level is normal because of increased globulin levels. These changes, however, can be detected if one measures the *albumin/globulin ratio*. Normally this ratio exceeds 1.0. The diseases just described that selectively affect albumin levels are associated with lesser ratios. Increased total protein levels, particularly the globulin fraction, occur with multiple myeloma and other gammopathies.

*Serum protein electrophoresis* (SPEP) can separate the various components of blood protein into bands or zones according to their electrical charge. Several well-established electrophoretic patterns have been identified and can be associated with specific diseases (Table 32).

**TABLE 32** Protein electrophoresis patterns in specific diseases

Pattern	Electrophoresis	Disease
Acute reaction	↓ Albumin ↑ Alpha <sub>2</sub> globulin	Acute infections, tissue necrosis, burns, surgery, stress, myocardial infarction
Chronic inflammatory	sl. ↓ Albumin sl. ↑ Gamma globulin N Alpha <sub>2</sub> globulin	Chronic infection, granulomatous diseases, cirrhosis, rheumatoid- collagen diseases
Nephrotic syndrome	↓↓ Albumin ↑↑ Alpha <sub>2</sub> globulin N ↑ Beta globulin	Nephrotic syndrome
Far-advanced cirrhosis	↓ Albumin ↑ Gamma globulin Incorporation of beta and gamma peaks	Far-advanced cirrhosis
Polyclonal gamma globulin elevation	↑↑ Gamma globulin with a broad peak	Cirrhosis, chronic infection, sarcoidosis, tuberculosis, endocarditis, rheumatoid- collagen diseases
Hypogamma- globulinemia	↓ Gamma globulin with normal other globulin levels	Light-chain multiple myeloma
Monoclonal gammopathy	Thin spikes in the beta (IgA, IgM) and gamma globulins	Myeloma, Waldenström macroglobulinemia, gammopathies

↓, Decreased; ↑, increased; sl. ↓, slightly decreased; sl. ↑, slightly increased; ↓↓, greatly decreased; ↑↑, greatly increased; N, normal.

If a spike is detected, *immunofixation electrophoresis (IFE)* can be done. In general, polyclonal spikes are associated with infectious or inflammatory diseases, whereas monoclonal-specific spikes are often neoplastic. IFE is used to indicate deficiencies or excesses as seen with macroglobulinemia, monoclonal



gammopathy of undetermined significance, and multiple myeloma. Specific monoclonal protein studies can be performed on the urine or the blood. Monoclonal immunoglobulin heavy chain (gamma, alpha, mu, delta, or epsilon) and/or light chains (kappa or lambda) can be identified.

Protein electrophoresis is also used to evaluate the major protein fractions found in urine. Normally only a small amount of albumin is seen. Urinary protein electrophoresis is useful in classifying the type of renal damage, if present. These electrophoresis techniques can be applied to the CSF or any other body fluid.

### Interfering factors

- Prolonged application of a tourniquet can increase both fractions of total proteins.
- Sampling of peripheral venous blood proximal to an IV administration site can result in an inaccurately low protein level. Likewise, massive IV infusion of crystalloid fluid can result in acute hypoproteinemia.
- ✚ Drugs that may cause *increased* protein levels include anabolic steroids, androgens, corticosteroids, dextran, growth hormone, insulin, phenazopyridine, and progesterone.
- ✚ Drugs that may cause *decreased* protein levels include ammonium ions, estrogens, hepatotoxic drugs, and oral contraceptives.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting or special preparation is required.

#### During

##### Blood

- Collect a venous blood sample in a gold-top tube.
- Indicate on the laboratory slip if the patient has received any vaccinations or immunizations within the past 6 months. Also, list any drugs that may affect test results.

##### Urine

- PT** Instruct the patient to begin the 24-hour collection after voiding. Follow guidelines on inside front cover.

#### After

- Apply pressure to the venipuncture site.

## Abnormal findings

▲ **Increased blood monoclonal immunoglobulins**

- Multiple myeloma
- Waldenström macroglobulinemia

▲ **Increased blood polyclonal immunoglobulins**

- Amyloidosis
- Autoimmune diseases
- Chronic infection or inflammation
- Chronic liver disease

▲ **Increased urine monoclonal immunoglobulins**

- Multiple myeloma
- Waldenström macroglobulinemia

See also [Table 32](#).

## notes

## protein C, protein S

**Type of test** Blood

### Normal findings

Protein S: 60%-130% of normal activity

Protein C: 70%-150% of normal activity

Protein C levels are lower in females and decrease with age in both males and females.

### Test explanation and related physiology

The plasma coagulation system is tightly regulated between thrombosis and fibrinolysis. This precise regulation is important. The protein C–protein S system is an important regulator of coagulation. Protein C inhibits the regulation of activated factor VIII and factor V (see [Figure 10](#), p. 251). This function of protein C is enhanced by protein S. Congenital deficiencies of these vitamin K–dependent proteins may cause spontaneous intravascular thrombosis. Furthermore, dysfunctional forms of the proteins result in a hypercoagulable state. In addition, nearly 50% of hypercoagulable states are caused by the presence of a factor V (factor V Leiden, p. 397) that is resistant to protein C inhibition. Acquired deficiencies are less commonly symptomatic.

These proteins are vitamin K dependent and are decreased in patients who are taking warfarin (Coumadin), in liver diseases, and in severe malnutrition. Because complement regulatory proteins are acute phase reactants, autoimmune diseases and other inflammatory diseases are associated with increased binding of protein S causing an acquired protein S deficiency. Affected patients may experience hypercoagulable events. Measurement of plasma-free protein S antigen is performed as the initial testing for protein S deficiency.

### Interfering factors

- Decreased protein C may occur in the postoperative state.
- Pregnancy or the use of exogenous sex hormones is associated with decreases in proteins C and S.
- Active clotting states, such as vein thrombosis, can lower levels of proteins S and C.
- ✚ Drugs that can *decrease* levels include vitamin K inhibitors (e.g., Coumadin).

## 752 protein C, protein S

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: no
  - Blood tube commonly used: blue
  - If more than one blood test is to be obtained, draw the blood for protein C or S second to avoid contamination with tissue thromboplastin that may occur in the first tube. If only blood for protein C or S is being drawn, draw a red-top tube first (and throw it away), and then draw the blood for this study in a blue-top tube (two-tube method of blood draw).
- PT** If the patient is found to be deficient in either protein, encourage the patient's family members to be tested, because they may be similarly affected.

### Abnormal findings

#### ▼ Decreased levels

Congenital deficiency of protein C or protein S  
DIC  
Hypercoagulable states  
Pulmonary emboli  
Arterial or venous thrombosis  
Vitamin K deficiency because of drugs or malnutrition  
Malignancy  
Autoimmune diseases  
Inflammation  
Warfarin

#### notes

## prothrombin time (PT, Protime, International normalized ratio [INR])

**Type of test** Blood

### Normal findings

(These depend on reagents used for PT.)

11.0-12.5 seconds; 85%-100%

Full anticoagulant therapy: > 1.5-2 times control value; 20%-30%

INR: 0.8-1.1

### Possible critical values

20 seconds

INR: > 5.5

### Test explanation and related physiology

The PT is used to evaluate the adequacy of the extrinsic system and common pathway in the clotting mechanism. The PT measures the clotting ability of factors I (fibrinogen), II (prothrombin), V, VII, and X. When these clotting factors exist in deficient quantities, PT is prolonged. Many diseases and drugs are associated with decreased levels of these factors. These include the following:

- *Hepatocellular liver disease* (e.g., cirrhosis, hepatitis, and neoplastic invasive processes). Factors I, II, V, VII, IX, and X are produced in the liver. With severe hepatocellular dysfunction, synthesis of these factors will not occur.
- *Obstructive biliary disease* (e.g., bile duct obstruction secondary to tumor or gallstones or intrahepatic cholestasis secondary to sepsis or drugs). As a result of the biliary obstruction, the bile necessary for fat absorption fails to enter the gut, and fat malabsorption results. Vitamins A, D, E, and K are fat soluble and also are not absorbed. Because the synthesis of factors II, VII, IX, and X depends on vitamin K, these factors will not be adequately produced and serum concentrations will fall.
- *Parenchymal (hepatocellular) liver disease* can be differentiated from obstructive biliary disease by determination of the patient's response to parenteral vitamin K administration. If PT returns to normal after 1 to 3 days of vitamin K administration (10 mg intramuscularly twice a day), one can safely assume that the patient has obstructive biliary disease that is causing vitamin K malabsorption. If, on the other hand, PT does not return to normal with the vitamin K injections,

## 754 prothrombin time

one can assume that severe hepatocellular disease exists and that the liver cells are incapable of synthesizing the clotting factors no matter how much vitamin K is available.

- *Oral anticoagulant administration.* The coumarin derivatives dicumarol and warfarin (Coumadin, Panwarfin) are used to prevent coagulation in patients with thromboembolic disease. These drugs interfere with the production of vitamin K–dependent clotting factors, which results in a prolongation of PT, as previously described. The adequacy of coumarin therapy can be monitored by following the patient’s PT. The ideal INR must be individualized for each patient (Table 33).

PT test results are usually given in seconds, along with a control value. The control value usually varies somewhat from day to day because the reagents used may vary. The patient’s PT should be approximately equal to the control value. Some laboratories report PT values as percentages of normal activity, because the patient’s results are compared with a curve representing normal clotting time. Normally, the patient’s PT is 85% to 100%.

To have uniform PT results for physicians in different parts of the country and the world, the World Health Organization has recommended that PT results include the use of the *INR* value. The reported INR results are independent of the reagents or methods used. Most hospitals are now reporting PT times in both absolute and INR numbers.

Such factors as weight, body mass index, age, diet, and concurrent medications are known to affect warfarin dose requirements during anticoagulation therapy. Warfarin interferes with the regeneration of reduced vitamin K from oxidized vitamin K in the *vitamin K oxidoreductase (VKOR)* complex. A recently identified gene for the major subunit of *VKOR*, called *VKORC1*,

**TABLE 33 Preferred international normalized ratio (INR) according to indication for anticoagulation**

Indication	Preferred INR
Deep vein thrombosis prophylaxis	1.5-2.0
Orthopedic surgery	2.0-3.0
Deep vein thrombosis	2.0-3.0
Atrial fibrillation	2.0-3.0
Pulmonary embolism	2.5-3.5
Prosthetic valve prophylaxis	3.0-4.0

has been identified and may explain up to 44% of the variance in warfarin dose requirements. Furthermore, warfarin is metabolized in part by the cytochrome P-450 enzyme CYP2C9. The *CYP2C9\*2* and *CYP2C9\*3* genetic mutations have been shown to decrease the enzyme activity of these metabolizing enzymes, which has led to warfarin sensitivity and, in serious cases, bleeding complications. A *warfarin pharmacogenomic test panel* is available that can identify any mutations in the *VKORC1-1639*, *CYP2C9\*2*, or *CYP2C9\*3* genes. The warfarin pharmacogenomic test can be used as part of an algorithm to determine the best initial warfarin dose and does not replace the need for routine PT testing for the calculation of the INR.

*Point-of-care home testing* is now available for patients who require long-term anticoagulation with warfarin. Like glucose monitoring, a finger stick is performed. A drop of blood is placed on the testing strip and inserted into the handheld testing device. The PT and INR are provided in a few minutes. The treating physician is notified by phone and any therapeutic changes can be instigated the same day.

### Interfering factors

- Alcohol intake can increase PT levels.
- A high-fat diet may decrease PT levels.
- ✚ Drugs that may cause *increased* levels include allopurinol, aminosalicic acid, barbiturates, beta-lactam antibiotics, cephalosporins, cholestyramine, chloral hydrate, chlorpromazine, cimetidine, clofibrate, colestipol, ethyl alcohol, glucagon, heparin, methyldopa, neomycin, oral anticoagulants, propylthiouracil, quinidine, quinine, salicylates, and sulfonamides.
- ✚ Drugs that may cause *decreased* levels include anabolic steroids, barbiturates, chloral hydrate, digitalis, diphenhydramine, estrogens, griseofulvin, oral contraceptives, and vitamin K.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: light blue
- Obtain the blood specimen before the patient is given the daily dose of warfarin.
- Remember, hemostasis will be delayed if the patient is taking warfarin or if the patient has any coagulopathies.
- Patients taking warfarin will have their doses regulated by PT and INR values.

## 756 prothrombin time

- If the PT is greatly prolonged, evaluate the patient for bleeding tendencies.
- PT** Teach patients on warfarin to check themselves for bleeding.
- The anticoagulant effect of warfarin can be reversed by the slow parenteral administration of vitamin K.
- PT** Instruct patients on warfarin therapy not to take any other medications unless approved by their physician.

### Abnormal findings

#### ▲ Increased levels

Cirrhosis  
Hepatitis  
Vitamin K deficiency  
Salicylate intoxication  
Bile duct obstruction  
Coumarin ingestion  
DIC  
Massive blood transfusion  
Hereditary factor deficiency

#### notes



## pulmonary angiography (Pulmonary arteriography)

**Type of test** X-ray with contrast dye

### Normal findings

Normal pulmonary vasculature

### Test explanation and related physiology

Through an injection of a radiographic contrast material into the pulmonary arteries, pulmonary angiography permits visualization of the pulmonary vasculature. Angiography is used to detect pulmonary embolism when the lung scan yields inconclusive results. This study has mostly been replaced by CT of the chest (see p. 281).

*Bronchial angiography* can be done to identify bleeding sites in the lungs. For this procedure, catheters are placed transarterially into the orifice of bronchial arteries. Radiopaque material is then injected, and the arteries are visualized. If a bleeding site is identified, the site can be injected with a sclerosing agent to prevent further bleeding.

### Contraindications

- Patients with allergies to iodinated dye
- Patients who are pregnant, unless benefits outweigh risks
- Patients with bleeding disorders

### Potential complications

- Allergic reaction to iodinated dye
- Hypoglycemia or acidosis may occur in patients who are taking metformin and receive iodine dye.
- Cardiac arrhythmia
- Premature ventricular contractions during right-sided heart catheterization may lead to ventricular tachycardia and ventricular fibrillation.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- Ensure that written and informed consent for this procedure is obtained.
- PT** Inform the patient that a warm flush will be felt when the dye is injected.
- Check the patient for allergies to iodinated dyes.

## 758 pulmonary angiography

- Determine whether the patient has ventricular arrhythmias.
- Keep the patient NPO after midnight on the day of the test.
- Administer preprocedural medications as ordered. Atropine may be given to decrease secretions. Meperidine may be used for sedation and relaxation.

### During

- Note the following procedural steps:
  1. The patient is placed on an x-ray table in the supine position.
  2. A catheter is placed into the femoral vein and passed into the inferior vena cava.
  3. With fluoroscopic visualization, the catheter is advanced into the main pulmonary artery, where the dye is injected.
  4. X-ray images of the chest are immediately taken in timed sequence. This allows all vessels visualized by the injection to be photographed. If filling defects are seen in the contrast-filled vessels, pulmonary emboli are present.
  5. If *bronchial angiography* is performed, the femoral artery is cannulated instead of the vein.
- Note that this test is performed by a physician in approximately 1 hour.

**PT** During injection of dye, inform the patient that he or she will feel a burning sensation and flush throughout the body.

### After

- Observe the catheter insertion site for inflammation, hemorrhage, and hematoma.
- Assess the patient's vital signs for evidence of bleeding.
- Apply cold compresses to the puncture site if needed to reduce swelling or discomfort.

### Abnormal findings

Pulmonary embolism

Congenital and acquired lesions of the pulmonary vessels (e.g., pulmonary hypertension)

Tumor

notes

## pulmonary function tests (PFTs)

**Type of test** Airflow assessment

### Normal findings

Vary with the patient's age, sex, height, and weight

### Test explanation and related physiology

PFTs are performed to detect abnormalities in respiratory function and to determine the extent of any pulmonary abnormality. The main reasons for pulmonary function tests include the following:

- Preoperative evaluation of the lungs and pulmonary reserve
- Evaluation of the response to bronchodilator therapy
- Differentiation between restrictive and obstructive forms of chronic pulmonary disease. *Restrictive* defects occur when ventilation is disturbed by a limitation in chest expansion. Inspiration is primarily affected. *Obstructive* defects occur when ventilation is disturbed by an increase in airway resistance. Expiration is primarily affected.
- Determination of the diffusing capacity of the lungs ( $D_L$ ). Rates are based on the difference in concentration of gases in inspired and expired air.
- Performance of inhalation tests in patients with allergies

PFTs include spirometry, measurement of airflow rates, and calculation of lung volumes and capacities. Exercise pulmonary stress testing can also be performed to provide data concerning the patient's pulmonary reserve.

On the basis of age, height, weight, race, and sex, normal values for the volumes and flow rates can be predicted. If the actual values are greater than 80% of predicted values, the results are considered normal. Spirometry provides information about obstruction or restriction of airflow. If airflow rates are significantly diminished (< 60% of normal), spirometry can be repeated after bronchodilators are administered by nebulizer.

PFTs include determination of the following:


- *Forced vital capacity (FVC)* is the amount of air that can be forcefully expelled from a maximally inflated lung position. This volume is decreased below the expected value in obstructive and restrictive pulmonary diseases.
- *Forced expiratory volume in 1 second ( $FEV_1$ )* is the volume of air expelled during the first second of the FVC. In patients with obstructive disease, airways are narrowed, and resistance to flow is high. Therefore not as much air can be expelled in

## 760 pulmonary function tests

$FEV_1$  will be reduced below the predicted value. In restrictive lung disease,  $FEV_1$  is decreased not because of airway resistance but because the amount of air originally inhaled is less. One should therefore measure the  $FEV_1/FVC$  ratio. A normal value of 80% is found in patients with restrictive lung disease. In obstructive lung disease, this ratio is considerably less than 80%.

- *Maximal midexpiratory flow (MMEF)* is the maximal rate of airflow through the pulmonary tree during forced expiration. This is also called *forced midexpiratory flow*. MMEF is reduced below expected values in obstructive diseases and normal in restrictive diseases.
- *Maximal volume ventilation (MVV)*, formerly called *maximal breathing capacity*, is the maximal volume of air that the patient can breathe in and out for 1 minute. MVV is decreased below the expected value in both restrictive and obstructive pulmonary disease.

A *comprehensive pulmonary function study* also may include evaluation of the following lung volumes and lung capacities, many of which are illustrated in [Figure 35](#).

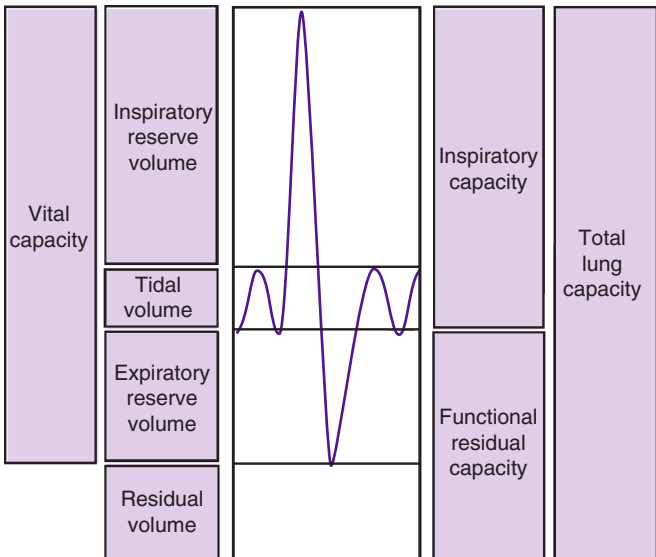

**FIG. 35** Relationship of lung volumes and capacities.

- *Tidal volume* ( $TV$  or  $V_T$ ) is the volume of air inspired and expired with each normal respiration.
- *Inspiratory reserve volume* ( $IRV$ ) is the maximal volume of air that can be inspired from the end of a normal inspiration. It represents forced inspiration over and beyond the tidal volume.
- *Expiratory reserve volume* ( $ERV$ ) is the maximal volume of air that can be exhaled after a normal expiration.
- *Residual volume* ( $RV$ ) is the volume of air remaining in the lungs after forced expiration.
- *Inspiratory capacity* ( $IC$ ) is the maximal amount of air that can be inspired after a normal expiration ( $IC = TV + IRV$ ).
- *Functional residual capacity* ( $FRC$ ) is the amount of air left in the lungs after a normal expiration ( $FRC = ERV + RV$ ).
- *Vital capacity* ( $VC$ ) is the maximal amount of air that can be expired after a maximal inspiration ( $VC = TV + IRV + ERV$ ).
- *Total lung capacity* ( $TLC$ ) is the volume to which the lungs can be expanded with the greatest inspiratory effort ( $TLC = TV + IRV + ERV + RV$ ).
- *Minute volume* ( $MV$ ), sometimes called *minute ventilation*, is the volume of air inhaled and exhaled per minute.
- *Dead space*. Dead space is the part of the tidal volume that does not participate in alveolar gas exchange. This would include the air in the trachea.
- *Forced expiratory flow*<sub>200-1200</sub> ( $FEF_{200-1200}$ ) is the airflow rate of expired air between 200 and 1200 mL during the FVC. This is the portion of the airflow curve that is the most affected by airway obstruction.
- *Forced expiratory flow*<sub>25-75</sub> ( $FEF_{25-75}$ ) is the airflow rate of expired air between 25% and 75% of the flow during the FVC. This is the part of the airflow curve that is the most affected by airway obstruction.
- *Peak inspiratory flow rate* ( $PIFR$ ) is the flow rate of inspired air during maximum inspiration. This is used to indicate large (trachea and bronchi) airway disease.
- *Peak expiratory flow rate* ( $PEFR$ ) is the maximum airflow rate during forced expiration.

Spirometry is the standard method for measuring most relative lung volumes; however, it is incapable of providing information about absolute volumes of air in the lung. Thus a different approach is required to measure  $RV$ ,  $FRC$ , and  $TLC$ . Two of the most common methods of obtaining information about these

## 762 pulmonary function tests

volumes are *body plethysmography* and *gas dilution tests*. Gas dilution or gas exchange studies are used to test *volume of isoflow* (*Viso V*), which is helpful in identifying early obstructive changes.

### Contraindications

- Patients who are in pain because of the inability for deep inspiration and expiration
- Patients who are unable to cooperate

### Procedure and patient care

#### Before

- PT** Explain the test to the patient.
- PT** Inform the patient that cooperation is necessary.
- PT** Instruct the patient not to use bronchodilators (if requested by health-care provider) or smoke for 6 hours before this test.
- PT** Tell the patient that the use of small-dose meter inhalers and aerosol therapy may be withheld before this study. Verify with health-care provider.
- Measure and record the patient's height and weight before this study to determine the predicted values.

#### During

- Note the following procedural steps:

#### *Spirometry and airflow rates*

1. The patient is taken to the pulmonary function laboratory.
2. The patient breathes through a sterile mouthpiece and into a spirometer to measure and record the values.
3. The patient is asked to inhale as deeply as possible and then forcibly exhale as much air as possible. This is repeated several times (usually two to three times). The two best are used for calculations.
4. From this, the machine computes FVC, FEV<sub>1</sub>, FEV<sub>1</sub>/FVC, PIFR, PEF, and MMEF.
5. The patient is asked to breathe in and out as deeply and frequently as possible for 15 seconds. The total volume breathed is recorded and multiplied by 4 to obtain the MVV.
6. The patient is asked to breathe in and out normally into the spirometer and then exhale forcibly from the end tidal volume expiration point. This provides the ERV.
7. The patient is asked to breathe in and out normally into the spirometer and then inhale forcibly from the end-tidal volume expiration point. This provides the IC.
8. The patient is asked to breathe in and out maximally. This is a measure of VC and the calculated TLC.

*Gas exchange/diffusing capacity of the lung*

- The  $D_L$  of CO is usually measured by having the patient inhale a CO mixture.
- $D_L$ CO is calculated with an analysis of the amount of CO exhaled compared with the amount inhaled.

*Inhalation tests (bronchial provocation studies)*

- These tests also may be performed during pulmonary function studies to establish a cause-and-effect relationship in some patients with inhalant allergies.
- The *Provocholine challenge* test is typically used to detect the presence of hyperactive airway diseases. This test would not be indicated for a patient with asthma.
- Care is taken during this challenge test to reverse any severe bronchospasm with prompt administration of an inhalant bronchodilator.

**After**

- Note that patients with severe respiratory problems are occasionally exhausted after the testing and will need rest.

**Abnormal findings**

Pulmonary fibrosis	Inhalation pneumonitis
Interstitial lung diseases	Postpneumonectomy
Tumor	Bronchiectasis
Chest wall trauma	Airway infection
Emphysema	Pneumonia
Chronic bronchitis	Neuromuscular disease
Asthma	Hypersensitivity bronchospasm

**notes**

**pyelography** (Intravenous pyelography [IVP], Excretory urography [EUG], Intravenous urography [IUG, IVU], Retrograde pyelography, Antegrade pyelography)

**Type of test** X-ray with contrast dye

### Normal findings

Normal size, shape, and position of the kidneys, renal pelvis, ureters, and bladder

Normal kidney excretory function as evidenced by the length of time for passage of contrast material through the kidneys

### Test explanation and related physiology

Pyelography is an x-ray study that uses radiopaque contrast material to visualize the kidneys, renal pelvis, ureters, and bladder. The contrast can be injected intravenously, through a catheter placed into the ureter (retrograde), or through a catheter placed into the proximal renal collecting system (antegrade).

IVP is indicated for patients with

- Pain compatible with urinary stones
- Blood in the urine
- Proposed pelvic surgery to locate the ureters
- Trauma to the urinary system
- Urinary outlet obstruction
- A suspected kidney tumor

For IVP, dye is injected intravenously, filtered out at the kidney by the glomeruli, and then passed through the renal tubules. X-ray images taken at set intervals over the next 30 minutes will show passage of the dye material through the kidneys and ureters and into the bladder.

If the artery leading to one of the kidneys is blocked, the dye cannot enter that part of the renal system, and that kidney or part thereof will not be visualized. If the artery is partially blocked, the length of time required for the appearance of the contrast material will be prolonged.

With primary glomerular disease (e.g., glomerulonephritis), the glomerular filtrate is reduced, which causes a reduction in the quantity of dye filtered. Therefore it requires more time for enough dye to enter the kidney filtrate and allow for renal opacification. As a result, kidney visualization is delayed. This indicates an estimate of renal function.

Defects in dye filling of the kidney can indicate renal tumors or cysts. Often intrinsic tumors, stones, extrinsic tumors, and



scarring can partially or completely obstruct the flow of dye through the collecting system (pelvis, ureters, bladder). CT scanning (see p. 272), however, is the method of choice for diagnosing urolithiasis.

If the obstruction has been of sufficient duration, the collecting system proximal to the obstruction will be dilated (hydronephrosis). Retroperitoneal and pelvic tumors, aneurysms, and enlarged lymph nodes also can produce extrinsic compression and distortions of the opacified collecting system.

IVP is used to assess a patient for congenital absence or malposition of the kidneys. Horseshoe kidneys (connection of the two kidneys), double ureters, and pelvic kidneys are typical congenital abnormalities.

*Retrograde pyelography* refers to radiographic visualization of the urinary tract through ureteral catheterization and the injection of contrast material. The ureters are catheterized during cystoscopy. A radiopaque material is injected into the ureters, and x-ray images are taken. Retrograde pyelography is helpful in radiographically examining the ureters when a ureter is obstructed.

*Antegrade pyelography* provides visualization of the renal pelvis for accurate placement of nephrostomy tubes. This study is used to identify the upper collecting system in an obstructed kidney to be used as a map for accurate percutaneous placement of a nephrostomy tube. This is performed on patients who have an obstruction of the ureter and hydronephrosis.

### Contraindications

- Patients who are allergic to iodinated dyes
- Patients who are severely dehydrated, because this can cause renal shutdown and failure
- Patients with renal insufficiency, as evidenced by a blood urea nitrogen value greater than 40 mg/dL because the iodinated nephrotoxic dye can worsen kidney function
- Patients with multiple myeloma because the iodinated nephrotoxic dye can worsen renal function
- Patients who are pregnant unless the benefits outweigh the risks of radiation exposure to the fetus

### Potential complications

- Allergy to iodine dye
- Infiltration of contrast dye
- Renal failure. This occurs most often in elderly patients who are chronically dehydrated before the dye injection.

## 766 pyelography

- Hypoglycemia or acidosis may occur in patients who are taking metformin (Glucophage) and receive iodine dye.
- Hemorrhage at the needle puncture site during antegrade pyelography, because the kidney is highly vascular
- Complications associated with *retrograde pyelography* include the following:
  1. Urinary tract infection
  2. Sepsis by seeding the bloodstream with bacteria from infected urine
  3. Perforation of the bladder or ureter
  4. Hematuria
  5. Temporary obstruction to ureter caused by ureteral edema

### Interfering factors

- Fecal material, gas, or barium in the bowel may obscure visualization of the renal system.
- Abnormal renal function studies may prevent adequate visualization of the urinary tract.
- Retained barium from previous studies may obscure visualization. Studies using barium (e.g., barium enema) should be scheduled after an IVP.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. Inform the patient that several x-ray images will be taken over 30 minutes. See p. xviii for radiation exposure and risks.
- Obtain informed consent if required by the institution.
- Check the patient for allergies to iodinated dye and shellfish.
- Give the patient a laxative or a cathartic, if ordered, the evening before the test.
- PT** Inform the patient of the required food and fluid restrictions. Some institutions prefer abstinence from solid foods for 8 hours before testing. Some allow a clear-liquid breakfast on the test day.
- Ensure adequate hydration for the patient (IV or oral) before and after the test to avoid dye-induced renal failure.
- Note that pediatric patients will have decreased fasting times, as ordered on an individual basis.
- Note that elderly and debilitated patients should have fasting times indicated specifically for them.
- Note that patients receiving high rates of IV fluids may have infusion rates decreased for several hours before the study to increase the concentration of the dye within the urinary system.

- Assess the patient's blood urea nitrogen and creatinine levels. Renal function could deteriorate as a result of the dye injection.
- If the antegrade or retrograde pyelography will be performed with the patient under general anesthesia, follow routine general anesthesia precautions. Keep the patient NPO after midnight on the day of the test. Fluids may be given intravenously.

### During

- Note the following procedural steps for *IVP*:
  1. The patient is taken to the radiology department and placed in the supine position.
  2. A plain image of the abdomen (KUB) is taken to ensure that no residual stool obscures visualization of the renal system. This also screens for calculi in the renal collecting system.
  3. Skin testing for iodine allergy is often done.
  4. A peripheral IV line is started (if not in place), and a contrast dye (e.g., Hypaque, Renografin) is given.
  5. X-ray images are taken at specific times, usually at 1, 5, 10, 15, 20, and 30 minutes and sometimes longer, to follow the course of the dye from the cortex of the kidney to the bladder.
  6. The patient is taken to the bathroom and asked to void.
  7. A postvoiding image is taken to visualize the empty bladder.
- PT** Inform the patient that the dye injection often causes a transitory flushing of the face, a feeling of warmth, a salty taste in the mouth, or even transient nausea. Initial IV needle placement and lying on a hard x-ray table are the only other discomforts associated with IVP.
- Note the following procedural steps for *retrograde pyelography*:
  1. The ureteral catheters are passed into the ureters by means of cystoscopy (see p. 309).
  2. Radiopaque contrast material is injected into the ureteral catheters, and x-ray images are taken.
- Note the following procedural steps for *antegrade pyelography*:
  1. The renal pelvis is localized by means of ultrasound.
  2. Under local anesthesia, a thin-walled needle is advanced into the lumen of the renal pelvis.
  3. Contrast material is injected and x-ray images are obtained.
  4. The nephrostomy tube is placed over guidewires and its position is affirmed by repeating the x-rays.

## 768 pyelography

**PT** Inform the patient that antegrade or retrograde pyelography is uncomfortable. If awake, the patient will feel pressure and an urge to void.

### After

**PT** Maintain adequate oral or IV hydration for several hours after the IVP to counteract fluid depletion caused by the test preparation. Encourage fluid intake.

**PT** Assess the patient's urinary output. A decreased output may be an indication of renal failure.

**PT** Evaluate elderly and debilitated patients for weakness because of the combination of fasting and catharsis necessary for test preparation. Instruct these patients to ambulate only with assistance.

- Note the color of the urine; a pink tinge is typically present. Report bright red blood or clots to the physician.
- See p. xxi for appropriate interventions concerning care for patients with iodine allergy.

### Abnormal findings

Pyelonephritis

Glomerulonephritis

Kidney tumor

Renal hematoma

Renal laceration

Cyst or polycystic disease of the kidney

Congenital abnormality of the urologic tract

Renal or ureteral calculi

Trauma to the kidneys, ureters, or bladder

Tumor of the collecting system

Hydronephrosis

Extrinsic compression of the collecting system (e.g., caused by tumor, aneurysm)

Bladder tumor

Prostate enlargement (male)

### notes

## rabies-neutralizing antibody test

**Type of test** Blood

### Normal findings

< 1:5

### Test explanation and related physiology

Identification and documentation of the presence of rabies virus-neutralizing antibody are important for veterinary health-care workers and others who may be or may have been exposed to the rabies virus. This test may be performed on patients who are at great risk for animal bites and have received the human diploid cell rabies vaccine (HDCV). A rabies titer of greater than 1:5 is considered to be protective.

The rabies virus antibody test is also used in making the diagnosis of rabies in a patient suspected of having been exposed to the virus. A fourfold rise in antibody titer over several weeks in a person not previously exposed to the HDCV indicates rabies exposure.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

Exposure to rabies vaccine

Recent bite exposure to rabies virus

Active rabies in a patient or animal

notes

**red blood cell count** (RBC count, Erythrocyte count)**Type of test** Blood**Normal findings**(RBC  $\times 10^6/\mu\text{L}$  or RBC  $\times 10^{12}/\text{L}$  [SI units])

Adult/elderly

Male: 4.7-6.1

Female: 4.2-5.4

Children

Newborn: 4.8-7.1

2-8 weeks: 4-6

2-6 months: 3.5-5.5

6 months-1 year: 3.5-5.2

1-6 years: 4-5.5

6-18 years: 4-5.5

**Test explanation and related physiology**

This test is a count of the number of circulating RBCs in 1  $\text{mm}^3$  of peripheral venous blood. The RBC count is routinely performed as part of a complete blood count (CBC).

Intravascular RBC trauma such as that caused by artificial heart valves or peripheral vascular atherosclerotic plaques shortens the life of the RBC. An enlarged spleen such as that caused by portal hypertension or leukemia may inappropriately destroy and remove normal RBCs from circulation.

Normal RBC values vary according to gender and age. Women tend to have lower values than men, and RBC counts tend to decrease with age. When the value is decreased lower than the range of the expected normal value, the patient is said to be anemic. Low RBC values are caused by decreased bone marrow production (e.g., myelofibrosis, leukemia, renal disease, or dietary deficiencies), increased blood loss (e.g., bleeding), or increased RBC destruction (hemolysis).

RBC counts greater than normal can be physiologically induced as a result of the body's requirements for greater oxygen-carrying capacity (e.g., at high altitudes). Diseases that produce chronic hypoxia (e.g., congenital heart disease) also provoke this physiologic increase in RBCs. Polycythemia vera is a neoplastic condition causing uncontrolled production of RBCs.

## Interfering factors

- Normal decreases are seen during pregnancy because of normal body fluid increases and dilution of the RBCs.
  - Persons living at high altitudes have increased RBCs.
  - Hydration status: Dehydration factitiously increases the RBC count, and overhydration decreases the RBC count.
- PT** Drugs that may cause *increased* RBC levels include erythropoietin and gentamicin.
- PT** Drugs that *decrease* RBC levels are many, including those that decrease marrow production or those that cause hemolysis.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender
- Thoroughly mix the blood with the anticoagulant by tilting the tube.
- Avoid hemolysis.

## Abnormal findings

### ▲ Increased levels

High altitude  
 Congenital heart disease  
 Polycythemia vera  
 Dehydration or hemoconcentration  
 Cor pulmonale  
 Pulmonary fibrosis  
 Thalassemia trait  
 Severe chronic obstructive pulmonary disease

### ▼ Decreased levels

Hemorrhage  
 Hemolysis  
 Anemia  
 Hemoglobinopathy  
 Advanced cancer  
 Bone marrow failure  
 Leukemia  
 Lymphoma  
 Antineoplastic chemotherapy  
 Chronic illness  
 Renal disease  
 Overhydration  
 Multiple myeloma  
 Pernicious anemia  
 Rheumatoid disease  
 Subacute endocarditis  
 Pregnancy  
 Dietary deficiency  
 Prosthetic valves

**red blood cell indices** (RBC indices, Blood indices)**Type of test** Blood**Normal findings**

Mean corpuscular volume (MCV)

Adult/elderly/child: 80-95 fL

Newborn: 96-108 fL

Mean corpuscular hemoglobin (MCH)

Adult/elderly/child: 27-31 pg

Newborn: 32-34 pg

Mean corpuscular hemoglobin concentration (MCHC)

Adult/elderly/child: 32-36 g/dL (or 32%-36%)

Newborn: 32-33 g/dL (or 32%-33%)

Red blood cell distribution width (RDW)

Adult: 11%-14.5%

**Test explanation and related physiology**

The RBC indices provide information about the size (MCV and RDW), hemoglobin content (MCH), and hemoglobin concentration (MCHC) of RBCs. This test is routinely performed as part of a CBC. The results of the RBC, hematocrit, and hemoglobin tests are necessary to calculate the RBC indices. When investigating anemia, it is helpful to categorize the anemia according to the RBC indices, as shown in [Table 34](#). Cell size is indicated by the terms *normocytic*, *microcytic*, and *macrocytic*. Hemoglobin content is indicated by the terms *normochromic*, *hypochromic*, and *hyperchromic*. Additional information about RBC size, shape, color, and intracellular structure is described in the blood smear study (see p. 149).

**Mean corpuscular volume**

The MCV is a measure of the average volume, or size, of a single RBC and is therefore used in classifying anemias. When the MCV value is increased, the RBC is said to be abnormally large, or *macrocytic*. This is most frequently seen in megaloblastic anemias (e.g., vitamin B<sub>12</sub> or folic acid deficiency). When the MCV value is decreased, the RBC is said to be abnormally small, or *microcytic*. This is associated with iron-deficiency anemia or thalassemia.



**TABLE 34** Categorization of anemia according to red blood cell indices**Normocytic,<sup>a</sup> normochromic<sup>b</sup> anemia**

Iron deficiency (detected early)  
 Chronic illness (e.g., sepsis, tumor)  
 Acute blood loss  
 Aplastic anemia (e.g., total body therapeutic irradiation)  
 Acquired hemolytic anemias (e.g., from a prosthetic cardiac valve)  
 Renal disease (because of the loss of erythropoietin)

**Microcytic,<sup>c</sup> hypochromic<sup>d</sup> anemia**

Iron deficiency (detected late)  
 Thalassemia  
 Lead poisoning

**Microcytic, normochromic anemia**

Chronic illnesses

**Macrocytic,<sup>e</sup> normochromic anemia**

Vitamin B<sub>12</sub> or folic acid deficiency  
 Phenytoin ingestion  
 Chemotherapy  
 Some myelodysplastic syndromes  
 Myeloid leukemia  
 Ethanol toxicity  
 Thyroid dysfunction

<sup>a</sup> Normocytic—normal red blood cell (RBC) size.

<sup>b</sup> Normochromic—normal color (normal hemoglobin content).

<sup>c</sup> Microcytic—smaller than normal RBC size.

<sup>d</sup> Hypochromic—less than normal color (decreased hemoglobin content).

<sup>e</sup> Macrocytic—larger than normal RBC size.

**Mean corpuscular hemoglobin**

The MCH is a measure of the average amount of hemoglobin in an RBC. MCH is calculated as follows:

$$\text{MCH} = \frac{\text{Hemoglobin (g/dL)} \times 10}{\text{RBC (million/mm}^3\text{)}}$$

Because macrocytic cells generally have more hemoglobin and microcytic cells have less hemoglobin, the causes for these values closely resemble those for the MCV value.

## 774 red blood cell indices

### Mean corpuscular hemoglobin concentration

The MCHC is a measure of the average concentration or percentage of hemoglobin in a single RBC. MCHC is calculated as follows:

$$\text{MCHC} = \frac{\text{Hemoglobin (g/dL)} \times 100}{\text{Hematocrit (\%)}}$$

When values are decreased, the cell has a deficiency of hemoglobin and is said to be *hypochromic* (frequently seen in iron-deficiency anemia and thalassemia). When values are normal, the anemia is said to be *normochromic* (e.g., hemolytic anemia). RBCs cannot be considered *hyperchromic*. Only 37 g/dL of hemoglobin can fit into the RBC. Alterations in RBC shape (spherocytosis), RBC agglutination, and a hemolyzed specimen may cause automated counting machines to indicate MCHC levels higher than normal.

### Red blood cell distribution width

The RDW is an indication of the variation in RBC size. It is calculated with a machine by using the MCV and RBC values. Variations in the width of RBCs may be helpful when classifying certain types of anemia. The RDW is essentially an indicator of the degree of *anisocytosis*, a blood condition characterized by RBCs of variable and abnormal size.

### Interfering factors

- Abnormal RBC size may affect indices.
- Extremely elevated white blood cell counts may affect RBC indices.
- Large RBC precursors (e.g., reticulocytes; see p. 788) cause an abnormally high MCV.
- Marked elevation in lipid levels (> 2000 mg/dL) causes automated cell counters to indicate high hemoglobin levels. MCHC and MCH will be calculated falsely high.
- The presence of cold agglutinins also falsely elevates MCHC, MCH, and MCV.
- ✦ Drugs that may *increase* MCV results include azathioprine, phenytoin, and zidovudine.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

**Abnormal findings****▲ Increased MCV**

Liver disease  
 Antimetabolite therapy  
 Alcoholism  
 Pernicious anemia (vitamin  
 $B_{12}$  deficiency)  
 Folic acid deficiency

**▲ Increased MCH**

Macrocytic anemia

**▲ Increased MCHC**

Spherocytosis  
 Intravascular hemolysis  
 Cold agglutinins

**▲ Increased RDW**

Iron deficiency anemia  
 $B_{12}$  or folate deficiency anemia  
 Hemoglobinopathies (e.g., sickle cell disease)  
 Hemolytic anemias  
 Posthemorrhagic anemias

**▼ Decreased MCV**

Iron deficiency anemia  
 Thalassemia  
 Anemia caused by chronic  
 illness

**▼ Decreased MCH**

Microcytic anemia  
 Hypochromic anemia

**▼ Decreased MCHC**

Iron deficiency anemia  
 Thalassemia

notes

**renal biopsy** (Kidney biopsy)

**Type of test** Microscopic examination of tissue

**Normal findings**

No pathologic conditions

**Test explanation and related physiology**

Biopsy of the kidney affords microscopic examination of renal tissue. Renal biopsy is performed for the following purposes:

- To diagnose the cause of renal disease (e.g., poststreptococcal glomerulonephritis, Goodpasture syndrome)
- To detect primary and metastatic malignancy of the kidney
- To evaluate the degree of rejection that occurs after kidney transplantation

Renal biopsy is most often obtained percutaneously (Figure 36). During this procedure, a needle is inserted through the skin and into the kidney to obtain a sample of kidney tissue. The biopsy needle is more accurately placed when guided by computed tomography (CT) scanning, ultrasonography, or fluoroscopy.

**Contraindications**

- Patients with coagulation disorders because of the risk of excessive bleeding
- Patients with operable kidney tumors because tumor cells may be disseminated during the procedure
- Patients with hydronephrosis because the enlarged renal pelvis can be easily entered and cause a persistent urine leak requiring surgical repair
- Patients with urinary tract infections because the needle insertion may disseminate the active infection


**Potential complications**

- Hemorrhage from the highly vascular renal tissue
- Inadvertent puncture of the liver, lung, bowel, aorta, and inferior vena cava
- Infection when an open biopsy is performed

**Procedure and patient care****Before**

**PT** Explain the procedure to the patient.

- Ensure that written and informed consent for this procedure is obtained by the physician.

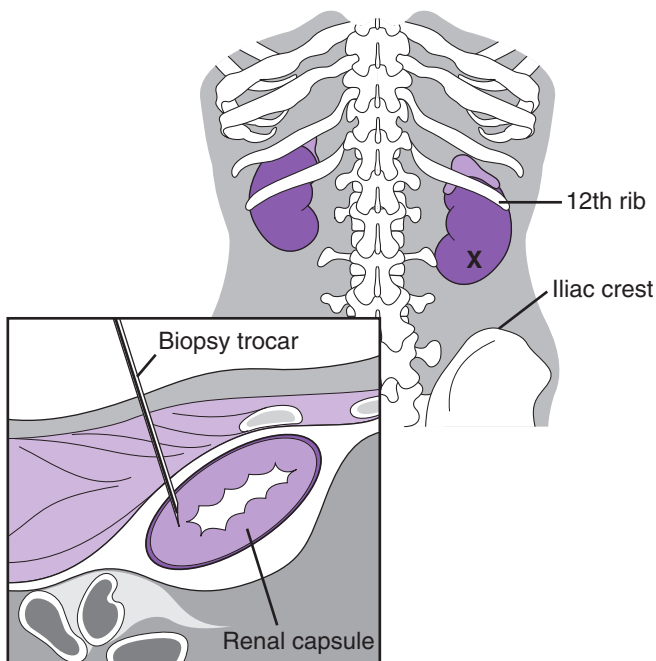

**FIG. 36** Renal biopsy. A biopsy needle is placed through the posterior lower chest wall and into the renal parenchyma, from which tissue is extracted.

- Keep the patient NPO (nothing by mouth) after midnight on the day of the test in the event that bleeding or inadvertent puncture of an abdominal organ may necessitate surgical intervention.
  - Assess the patient's coagulation studies.
  - Check the patient's hemoglobin and hematocrit values.
  - Note that the patient may need to be typed and crossmatched for blood in the event of severe hemorrhage.
- PT** Tell the patient that no sedative is required.
- Note that the needle stick may be done at the bedside.
  - If CT scanning or ultrasound guidance is used, the needle stick is done in the radiology or ultrasonography department.

## 778 renal biopsy

### During

- Note the following procedural steps:
  1. The patient is placed in a prone position with a sandbag or pillow under the abdomen to straighten the spine.
  2. Under sterile conditions, the skin overlying the kidneys is infiltrated with a local anesthetic (lidocaine).
  3. While the patient holds his or her breath to stop kidney motion, the physician inserts the biopsy needle into the kidney and takes a specimen.
  4. After this procedure is completed, the needle is removed, and pressure is applied to the site.
- Note that this procedure is performed by a physician in approximately 10 minutes.
- PT** Tell the patient that this procedure is uncomfortable but only minimally if enough lidocaine is used.

### After

- After the test, apply a pressure dressing.
- Turn the patient onto his or her back and keep him or her on bed rest for approximately 24 hours.
- Check the patient's vital signs, puncture site, and hematocrit values frequently during the 24-hour period.
- PT** Instruct the patient to avoid any activity that increases abdominal venous pressure (e.g., coughing).
- Assess the patient for signs and symptoms of hemorrhage (e.g., decrease in blood pressure, increase in pulse, pallor).
- Evaluate the patient's abdomen for signs of bowel or liver penetration (e.g., abdominal pain and tenderness, abdominal muscle guarding and rigidity, decreased bowel sounds).
- PT** Instruct the patient to avoid strenuous exercise or jolting of the kidney for at least 2 weeks.
- PT** Instruct the patient to report burning during urination or development of a fever. Either could indicate an infection.
- Inspect all urine specimens for gross hematuria. Usually the patient's urine will contain blood initially, but this generally will not continue after the first 24 hours.
- PT** Encourage the patient to drink large amounts of fluid to prevent clot formation and urine retention.

### Abnormal findings

Renal disease (e.g., poststreptococcal conditions, glomerulonephritis, lupus nephritis)

Primary and metastatic malignancy of the kidney

Rejection of kidney transplant

### notes

**renal scanning** (Kidney scan, Radiorenography, Radionuclide renal imaging, Nuclear imaging of the kidney, DSMA renal scan, DTPA renal scan, Captopril renal scan)

**Type of test** Nuclear scan

### **Normal findings**

Normal size, shape, and function of the kidney

### **Test explanation and related physiology**

Renal scans are used to indicate the perfusion, function, and structure of the kidneys. They are also used to indicate the presence of ureteral obstruction or renovascular hypertension. Because this study uses no iodinated dyes, it is safe to perform on patients who have iodine allergies or compromised renal function. Renal scans are used to monitor renal function in patients with known renal disease. This scan also plays a large part in the diagnosis of renal transplant rejection.

This nuclear procedure provides visualization of the urinary tract after intravenous (IV) administration of a radioisotope.

### **Renal blood flow (perfusion) scan**

This type of renal scan is used to evaluate the blood flow to each kidney. It is used to identify renal artery stenosis, renovascular hypertension, and rejection of renal transplant. Also, it is used to demonstrate hypervascular lesions (renal cell carcinoma) in the kidney.

Decreased gamma activity is noted in a kidney with arterial stenosis or renovascular hypertension. Decreased activity relative to the aorta is noted in a transplanted kidney that is experiencing rejection. Localized increased gamma activity is noted in a kidney that contains a hypervascular tumor (cancer).

### **Renal structural scan**

This type of renal scan is performed to outline the structure of the kidney to identify pathology that may alter normal anatomic structure (e.g., tumor, cyst, abscess). Congenital disorders (e.g., hypoplasia or aplasia of the kidney, malposition of the kidney) can also be detected. A filling defect in the renal parenchyma may indicate a tumor, cyst, abscess, or infarction. A horseshoe-shaped kidney, pelvic kidney, or absence of a kidney may be evident. Also, information concerning postrenal transplants can be obtained with this scan. Anatomic alterations in the parenchymal distribution of tracer may indicate transplant rejection.

### Renal function scan (Renogram)

Renal function can be determined by documenting the capability of the kidney to take up and excrete a particular radioisotope. A well-functioning kidney can be expected to rapidly assimilate the isotope and then excrete it. A poorly functioning kidney will not be able to take up the isotope rapidly or excrete it in a timely manner. Renal function can be monitored by serially repeating this test and comparing results. Each radioactive tracer is handled by the kidney in a different manner. Different renal functions can be tested according to which isotope is used:

- Technetium-99 m diethylenetriamine pentaacetic acid ( $^{99m}\text{Tc}$ -DTPA) measures glomerular filtration.
- Technetium-99m-mercaptoacetyltriglycine ( $^{99m}\text{Tc}$ -MAG3) measures both glomerular filtration and tubular cell secretion.

### Renal hypertension scan

This scan is used to identify the presence and location of renovascular hypertension. It usually uses angiotensin-converting enzyme (ACE) inhibitors such as captopril. The *captopril scan* (captopril renography/scintigraphy) determines the functional significance of a renal artery or arteriole stenosis. These scans may predict the response of the blood pressure to medical treatment, angioplasty, or surgery.

### Renal obstruction scan

This scan is performed to identify obstruction of the outflow tract of the kidney caused by obstruction of the renal pelvis, ureter, or bladder outlet. Ultrasound, CT scanning, or magnetic resonance imaging are preferable and more accurate for anatomic abnormalities, tumors, and cysts.

Often several of these scans are combined to obtain the maximum possible information about the renal system. A *triple renal study* may use all of these techniques to evaluate renal blood perfusion, structure, and excretion. Radionuclear scans are also helpful in the evaluation of arterial trauma.

### Contraindications

- Patients who are pregnant, unless the benefits outweigh the risks of fetal damage

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.



- Do not schedule a renal scan within 24 hours after an intravenous pyelogram.
- PT** Assure the patient that he or she will not be exposed to large amounts of radioactivity because only tracer doses of isotopes are used.
- Remind the patient to void before the scan.
- PT** Tell the patient that no sedation or fasting is required but that good hydration is essential.
- PT** Instruct the patient to drink two to three glasses of water before the scan.

### During

- Note the following procedural steps:
  1. The unседated, nonfasting patient is taken to the nuclear medicine department.
  2. A peripheral IV injection of radionuclide is given. It takes only minutes for the radioisotopes to be concentrated in the kidneys.
  3. While the patient assumes a supine or sitting position, a gamma ray scintigraphy camera is passed over the kidney area and records the radioactive uptake on film.
  4. Scans may be repeated at different intervals after the initial isotope injection.
  5. Unique features of the various scans:
    - a. For a *furosemide (Lasix) renal scan* or a *diuretic renal scan*, images are obtained for 10 to 20 minutes; then furosemide is administered through an IV, and another 20 minutes of images is obtained.
    - b. For the *captopril renal scan*, the patient is scanned after the administration of captopril.
    - c. For the *renal blood flow* and the *renal function scans*, scanning is started immediately after the injection.
    - d. For *structural renal scans*, the patient is asked to lie still for the entire time of the scan (30 minutes).
- Note that the duration of this test varies from 1 to 4 hours, depending on the specific information required. Perfusion scans are done in approximately 20 minutes and functional scans in less than 1 hour. Static structure scans require 20 minutes to 4 hours for completion.
- Note that this study is performed by a nuclear medicine technologist or physician.
- PT** Tell the patient that no pain or discomfort is associated with this procedure.
- PT** Inform the patient of the need to lie still during this study.

## 782 renal scanning

### After

- PT** Because only tracer doses of radioisotopes are used, inform the patient that no precautions need to be taken against radioactive exposure.
- PT** Tell the patient that the radioactive substance is usually excreted from the body within 6 to 24 hours. Encourage the patient to drink fluids.

### Abnormal findings

Urinary obstruction  
Pyelonephritis  
Renovascular hypertension  
Absence of kidney function  
Renal infarction  
Renal arterial atherosclerosis  
Glomerulonephritis  
Renal tumor  
Congenital abnormalities  
Renal trauma  
Transplant rejection  
Acute tubular necrosis  
Renal abscess  
Renal cyst

### notes

## renin assay, plasma (Plasma renin activity [PRA], Plasma renin concentration [PRC])

**Type of test** Blood

### Normal findings

#### Plasma renin assay

Adult/elderly

*Upright position, sodium depleted* (sodium-restricted diet)

Ages 20-39 years: 2.9-24 ng/mL/hr

> 40 years: 2.9-10.8 ng/mL/hr

*Upright position, sodium repleted* (normal sodium diet)

Ages 20-39 years: 0.6-4.3 ng/mL/hr

> 40 years: 0.6-3 ng/mL/hr

Child

0-3 years: < 16.6 ng/mL/hr

3-6 years: < 6.7 ng/mL/hr

6-9 years: < 4.4 ng/mL/hr

9-12 years: < 5.9 ng/mL/hr

12-15 years: < 4.2 ng/mL/hr

15-18 years: < 4.3 ng/mL/hr

#### Renal vein


Renin ratio of involved kidney to uninvolved kidney < 1.4

### Test explanation and related physiology

Renin is an enzyme released by the juxtaglomerular apparatus of the kidney into the renal veins in response to hyperkalemia, sodium depletion, decreased renal blood perfusion, or hypovolemia. Renin activates the renin-angiotensin system ([Figure 37](#)).

Renin is not actually measured in this test. Plasma renin activity (PRA) measures enzyme ability to convert angiotensinogen to angiotensin I and is limited by the availability of angiotensinogen. The PRA test actually measures, by radioimmunoassay, the rate of angiotensin I generation per unit time. This is a commonly used renin assay. The specimen must be drawn under ideal circumstances, handled by the local laboratory correctly, and transferred to the central laboratory in a timely manner. Even then, results may vary significantly.

The PRA test is a screening procedure for the detection of renal-based or renovascular hypertension. The PRA may be supplemented by other tests such as the renal vein renin assay. A determination of the PRA and a simultaneous measurement

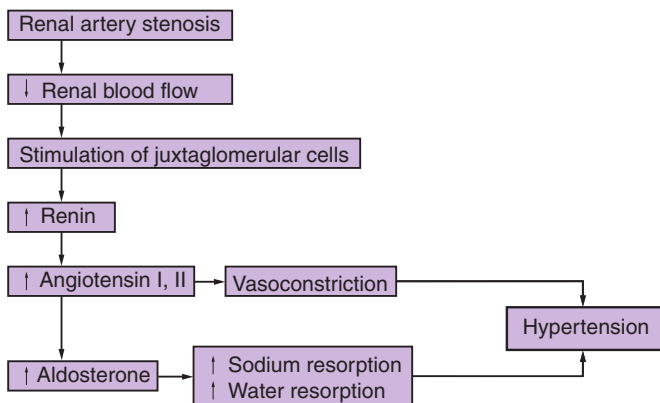

**FIG. 37** Physiology of renovascular hypertension.

of the plasma aldosterone (see p. 25) level are used in the differential diagnosis of primary versus secondary hyperaldosteronism. Patients with primary hyperaldosteronism (adrenal adenoma overproducing aldosterone or Conn syndrome) will have increased aldosterone production associated with suppressed renin activity. The aldosterone/renin ratio is greater than or equal to 20. Patients with secondary hyperaldosteronism (caused by renovascular occlusive disease or primary renal disease) will have increased levels of aldosterone and plasma renin.

*Renal vein assays* for renin are used to diagnose and lateralize renovascular hypertension—that is, hypertension that is related to inappropriately high renin levels from a diseased kidney or a hypoperfused kidney. The renal veins can be identified using injection of a radiopaque dye into the inferior vena cava. A catheter is placed into each renal vein, and blood is withdrawn from each vein. PRA is determined in each sample. If hypertension is caused by renal artery stenosis or renal pathology, the renal vein renin level of the affected kidney should be 1.5 or more times greater than that of the unaffected kidney or peripheral venous sample. If the levels are the same, the hypertension is not caused by a renovascular source.

The *renin stimulation test* can be performed to more clearly diagnose and distinguish primary and secondary hyperaldosteronism. In this test, PRA is obtained while the patient is in the recumbent position and on a low-salt diet. The PRA is then repeated with the patient on the same diet while the patient

is standing erect. In primary hyperaldosteronism, the blood volume is greatly expanded. A change in position or reduced salt intake will not result in decreased renal perfusion or sodium level. Therefore renin levels do not increase. In secondary hyperaldosteronism (or normal persons with essential hypertension), the renal perfusion decreases while in the upright position and sodium levels do decrease with decreased intake. Therefore renin levels increase.

The PRA is assessed as part of the *captopril test* (a screening test for renovascular hypertension). Patients with renovascular hypertension have greater falls in blood pressure and increases in PRA after administration of ACE inhibitors than do those with essential hypertension. For the captopril test, the patient receives an oral dose of captopril (ACE inhibitor) after a baseline PRA test, and blood pressure measurements are then taken. Subsequent blood pressure measurements and a repeat PRA test at 60 minutes are used for test interpretation. This is an excellent screening procedure to determine the need for a more invasive radiographic evaluation (e.g., digital subtraction renal arteriography or bilateral renal arteriography).

### Contraindications

- Patients who are allergic to iodinated dye

### Potential complications

- Allergic reaction to iodinated dye

### Interfering factors

- Renin levels are affected by pregnancy, salt intake, and licorice ingestion.
- Values are higher in patients on low salt diets and when the patient is in an upright position.
- Posture: Renin is increased in the erect position and decreased in the recumbent position.
- Values are higher early in the day. There is a diurnal variation in renin production.
- ✘ Spironolactone interferes with renin testing and should be discontinued 4 to 6 weeks before renin testing.
- ✘ Drugs that *increase* levels of renin include ACE inhibitors, antihypertensives, diuretics, estrogens, oral contraceptives, and vasodilators.
- ✘ Drugs that *decrease* renin levels include beta-blockers, clonidine, licorice, nonsteroidal antiinflammatory agents, potassium, and reserpine.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to maintain a normal diet with a restricted amount of sodium (~3 g/day) for 3 days before the test. A high-sodium diet causes a decrease in renin.
  - A urine sodium (see p. 837) may be obtained to normalize PRA to salt intake.
- PT** Instruct the patient to discontinue medications that may interfere with results for 2 to 4 weeks before the test as ordered by the physician.
  - Usually draw a fasting blood sample because renin values are higher in the morning.
- PT** For stimulation tests, instruct the patient to significantly reduce sodium intake for 3 days before testing.
  - Ensure that the patient stands or sits upright for 2 hours before the blood is drawn.
  - If a recumbent sample is ordered, have the patient remain in bed in the morning until the blood has been obtained.

### During

- The test is usually performed with the patient in an upright position.
- For the stimulation tests, the blood is drawn in the recumbent and upright positions.
- Release the tourniquet immediately before obtaining the blood specimen because stasis can lower renin levels.
- Collect a venous blood sample in a chilled lavender-top tube.
- Gently invert the blood tube to allow adequate mixing of the blood sample and the anticoagulant.
- Record the patient's position and dietary status and time of day on the laboratory slip.
- Place the tube of blood on ice and immediately send it to the laboratory.

### After

- Apply pressure to the venipuncture site.
- PT** Tell the patient that usually a normal diet may be resumed.

## Abnormal findings

### ▲ Increased levels

- Essential hypertension
- Malignant hypertension
- Renovascular hypertension
- Chronic renal failure
- Salt-losing gastrointestinal disease (e.g., diarrhea)
- Addison disease
- Renin-producing renal tumor
- Bartter syndrome
- Cirrhosis
- Hyperkalemia
- Hemorrhage
- Hypovolemia

### ▼ Decreased levels

- Primary hyperaldosteronism
- Steroid therapy
- Congenital adrenal hyperplasia
- Chronic renal impairment
- Hypervolemia
- Ectopic adrenocorticotrophic hormone syndrome

## notes

## reticulocyte count (Retic count)

**Type of test** Blood

### Normal findings

Reticulocyte count:

Adult/elderly/child: 0.5%-2%

Infant: 0.5%-3.1%

Newborn: 2.5%-6.5%

Reticulocyte index: 1.0

### Test explanation and related physiology

A reticulocyte is an immature RBC that can be readily identified under a microscope. Normally there are a small number of reticulocytes in the bloodstream. The reticulocyte count is a test for determining bone marrow function and evaluating erythropoietic activity. This test is also useful in classifying anemias.

Increased reticulocyte counts indicate that the marrow is putting an increased number of RBCs into the bloodstream, usually in response to anemia. A normal or low reticulocyte count in a patient with anemia indicates that the marrow response to the anemia by way of production of RBCs is inadequate and perhaps contributing to or the cause of the anemia (as in aplastic anemia, iron deficiency, vitamin B<sub>12</sub> deficiency, or depletion of iron stores). An elevated reticulocyte count found in patients with a normal hemogram indicates increased RBC production compensating for an ongoing loss of RBCs (hemolysis or hemorrhage).

To determine whether a reticulocyte count indicates an appropriate erythropoietic (RBC marrow) response in patients with anemia and a decreased hematocrit, one should calculate the *reticulocyte index*:

$$\text{Reticulocyte index} = \text{Reticulocyte count (in\%)} \times \frac{\text{Patient's hematocrit}}{\text{Normal hematocrit}}$$

The reticulocyte index in a patient with a good marrow response to the anemia should be 1.0. If it is lower than 1.0, even though the reticulocyte count is elevated, the bone marrow response is inadequate in its ability to compensate.

Determination of the *reticulocyte-specific hemoglobin content* (or *reticulocyte hemoglobin equivalent*) is a measure of the mean hemoglobin in reticulocytes. This test indicates the amount of



iron available for incorporation into hemoglobin over the previous 3 to 5 days. It is a very reliable test to identify iron deficiency, especially in children or in the face of complex other chronic diseases.

### Interfering factors

- Pregnancy may cause an increased reticulocyte count.
- RBCs containing Howell-Jolly bodies look like reticulocytes and can be miscounted by some automated counter machines to be reticulocytes and give a falsely high number of reticulocytes.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

### Abnormal findings

#### ▲ Increased levels

Hemolytic anemia  
 Sickle cell anemia  
 Hemorrhage (3-4 days later)  
 Postsplenectomy  
 Hemolytic disease of the newborn  
 Pregnancy  
 Leukemias  
 Recovery from nutritional anemias

#### ▼ Decreased levels

Pernicious anemia  
 Folic acid deficiency  
 Iron-deficiency anemia  
 Adrenocortical hypofunction  
 Aplastic anemia  
 Radiation therapy  
 Marrow failure  
 Anterior pituitary hypofunction  
 Chronic infection  
 Cirrhosis  
 Malignancy

notes

**rheumatoid factor (RF)****Type of test** Blood**Normal findings**

Negative (< 60 units/mL by nephelometric testing)  
(Elderly patients may have slightly increased values.)

**Test explanation and related physiology**

In rheumatoid arthritis (RA), abnormal immunoglobulin G (IgG) antibodies produced by lymphocytes in the synovial membranes act as antigens. Other IgG and IgM antibodies in the patient's serum react with the  $f_c$  component of the abnormal synovial antigenic IgG to produce immune complexes. These immune complexes activate the complement system and other inflammatory systems to cause joint damage. The reactive IgM is called *RF*. Tissues other than the joints, including those of the blood vessels, lungs, nerves, and heart, may also be involved in the autoimmune inflammation.

Tests for RF are directed toward identification of the IgM antibodies. Approximately 80% of patients with RA have positive RF titers. To be considered positive, RF must be found in a dilution of greater than 1:80; when RF is found in titers less than 1:80, such diseases as systemic lupus erythematosus, scleroderma, and other autoimmune conditions should be considered. Although the normal value is “no rheumatoid factor identifiable at low titers,” a small number of normal patients will have RF present at a very low titer. Furthermore, a negative RF does not exclude the diagnosis of RA. RF is not a useful disease marker because its presence does not disappear in patients who are experiencing a remission from the disease symptoms.

Other autoimmune diseases (see [Table 1](#), p. 87) such as systemic lupus erythematosus and Sjögren syndrome also may cause a positive RF test. RF is occasionally seen in patients with tuberculosis, chronic hepatitis, infectious mononucleosis, and subacute bacterial endocarditis as well.

**Interfering factors**

- Elderly patients often have false-positive results.
- Hemolysis or lipemia may be associated with false-positive results.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

## Abnormal findings

### ▲ Increased levels

RA

Other autoimmune diseases (e.g., systemic lupus erythematosus)

Chronic viral infection

Subacute bacterial endocarditis

Tuberculosis

Chronic hepatitis

Dermatomyositis

Scleroderma

Infectious mononucleosis

Leukemia

Cirrhosis

Syphilis

Renal disease

## notes

**ribosome P antibodies** (Ribosomal P Ab, Anti-ribosome P antibodies)**Type of test** Blood**Normal findings**

Negative

**Test explanation and related physiology**

This antibody test should not be confused with the antiextractable nuclear antibodies (antiribonucleoprotein antibody; see p. 75). Ribosome P antibodies are used as an adjunct in the evaluation of patients with lupus erythematosus (LE). Antibodies to ribosome P proteins are considered highly specific for LE and have been reported in patients with central nervous system (CNS) involvement (i.e., lupus psychosis). This antibody is therefore an aid in the differential diagnosis of neuropsychiatric symptoms in patients with LE. Because patients with LE may manifest signs and symptoms of CNS diseases, including neuropsychiatric symptoms, the presence of antibodies to ribosome P protein may be useful in the differential diagnosis of such patients. Most patients with LE do not have detectable levels of antibodies to ribosome P protein. But when they do, CNS involvement should be considered possible.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Check the venipuncture site for infection. Patients with autoimmune disease have compromised immune systems.

**Abnormal findings****▲ Increased levels**

Lupus erythematosus

**notes**

## rubella antibody test (German measles test)

**Type of test** Blood

### Normal findings

Method	Result	Interpretation
HAI titer	< 1:8	No immunity to rubella
HAI titer	> 1:20	Immune to rubella
LA	Negative	No immunity to rubella
ELISA IgM	< 0.9 IU/mL	No infection
ELISA IgM	> 1.1 IU/mL	Active infection
ELISA IgG	< 7 IU/mL	No immunity to rubella
ELISA IgG	> 10 IU/mL	Immune to rubella

*ELISA*, Enzyme-linked immunosorbent assay; *HAI*, hemagglutination inhibition; *LA*, latex agglutination.

### Possible critical values

Evidence of susceptibility in pregnant women with recent exposure to rubella

### Test explanation and related physiology

Screening for rubella antibodies is done to detect immunity to rubella. These tests detect the presence of IgG and/or IgM antibodies to rubella (the causative agent for German measles). They become elevated a few days to a few weeks after the onset of the rash, depending on what method of testing is used. IgM tends to disappear after about 6 weeks. IgG, however, persists at low but detectable levels for years. These antibodies become elevated in patients with active rubella infection or with past infections.

Children are vaccinated with rubella to prevent the effects of the disease and to minimize infection. Rubella testing documents immunity to rubella. Rubella immunity testing is suggested for all healthcare workers. Most importantly, however, it is done to verify the presence or absence of rubella immunity in pregnant women because congenital rubella infection in the first trimester of pregnancy is associated with congenital abnormalities of the fetus (heart defects, brain damage, deafness), abortion, or stillbirth. The term TORCH (toxoplasmosis, other, rubella,

## 794 rubella antibody test

cytomegalovirus, herpes) has been applied to infections with recognized detrimental effects on the fetus. All of these tests are discussed separately (see separate listings).

If the woman's titer is greater than 1:20, she is not susceptible to rubella. If the woman's titer is 1:8 or less, she has little or no immunity to rubella. A fourfold increase in HAI rubella titer from the acute to the convalescent titer indicates that the rash was caused by an active rubella infection. Alternatively, an IgM antibody titer could be done. If the titer is positive, recent infection has occurred. IgM titers appear 1 to 2 days after onset of the rash and disappear 5 to 6 weeks after infection.

Antirubella antibody testing is also used to diagnose rubella in infants (congenital rubella). IgM antirubella antibodies cannot pass through the placenta. If an infant has IgM antibodies, acute congenital or newborn rubella is suspected. Antibody testing is often used in children with congenital abnormalities that may have come from congenital rubella infection.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: no
  - Blood tube commonly used: red
- PT** Inform the patient when to return for a follow-up HAI titer if indicated.

### Abnormal findings

Active rubella infection

Previous rubella infection leading to immunity

### notes

## salivary gland nuclear imaging (Parotid gland nuclear imaging)

**Type of test** Nuclear medicine

### Normal findings

Normal function of the salivary gland; no tumor or duct obstruction

### Test explanation and related physiology

The ability of the epithelial cells of the salivary glands to transport large pertechnetate ions from the blood and to secrete them into the saliva provides the principle for imaging the salivary glands. The functional capabilities, structural integrity, and location of the glands can be assessed. Usually, the parotid gland alone is visualized. Occasionally, the submandibular glands can be seen.

Indications for salivary gland nuclear imaging include patients with the following:

- Xerostomia (dry mouth)
- Pain
- Tumors
- Possible parotid duct obstruction

By following the radionuclide immediately after injection, blood flow can be evaluated. In about 10 minutes after injection, gland function becomes obvious by uptake of the radionuclide into the gland. Five to 10 minutes later, one should see secretion of nuclear material into the mouth. Washout demonstrates complete salivary gland excretion. Usually the patient is asked to suck on a lemon to encourage rapid washout. This test can indicate inflammation, hypofunction, the location and character of tumors, and duct obstruction.

### Contraindications

- Patients who are pregnant unless the benefits outweigh the risks

### Interfering factors

- Rinsing the mouth before study may reduce excretion.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.

**PT** Tell the patient that no specific preparation is necessary.

- Make certain that the patient does not receive any thyroid-blocking agents within 48 hours of testing.

## 796 salivary gland nuclear imaging

### During

- Tc-99m pertechnetate is injected into the antecubital vein.
- Dynamic planar images are obtained immediately.
- Repeat images are obtained every 3 to 5 minutes for a total of 15 to 20 minutes.
- Three-dimensional images are often obtained by using single photon emission computed tomography (SPECT)/computed tomography (CT) imaging.
- Administer a salivary gland stimulant after completion of static images. Either lemon juice or a lemon slice should be swished in the mouth and then expectorated.
- Washout images are obtained 5 to 10 minutes after the salivary gland stimulant. The thyroid gland is included for reference or comparison.

### After

- PT** Assure the patient that the dose of radioactive technetium used in this test is minute and therefore harmless. No isolation or special urine precautions are needed.

### Abnormal findings

Sjögren syndrome

Benign mixed tumors or pleomorphic adenomas

Malignant lesions (e.g., adenocarcinomas, squamous cell carcinomas, undifferentiated and mixed carcinomas)

Salivary duct obstruction

### notes



## SARS viral testing

**Type of test** Blood; microscopic

### Normal findings

No SARS virus

### Test explanation and related physiology

Severe acute respiratory syndrome (SARS) is caused by a coronavirus (CoV). Routine testing for the SARS virus is not conducted unless a cluster of cases develops and health officials are able to rule out all other infectious agents.

The diagnosis can only be made with positive test results in the following situations with:

- One specimen tested on two occasions using the original clinical specimen on each occasion;
- Two clinical specimens from different sources (e.g., nasopharyngeal and stool); and
- Two clinical specimens collected from the same source on two different days (e.g., two nasopharyngeal aspirates).

Eight types of respiratory specimens may be collected for viral and/or bacterial diagnostics: (1) nasopharyngeal wash or aspirates, (2) nasopharyngeal swabs, (3) oropharyngeal swabs, (4) bronchoalveolar lavage, (5) tracheal aspirate, (6) pleural fluid tap, (7) sputum, and (8) postmortem tissue. A nasopharyngeal wash or aspirate is the specimen of choice for detection of most respiratory viruses.

Serum and blood (plasma) should be collected early in the illness. Both acute and convalescent serum specimens should be collected for antibody testing. To confirm or rule out SARS-CoV infection, it is important to collect convalescent serum specimens more than 28 days after the onset of illness.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Observe all universal precautions in handling the specimen.
- Observe strict isolation technique. This disease is contagious.

#### During

- To obtain a *nasopharyngeal wash/aspirate*, have the patient sit with the head tilted slightly backward. Instill 1 to 1.5 mL of nonbacteriostatic saline (pH 7.0) into one nostril. Flush a plastic catheter or tubing with 2 to 3 mL of saline. Insert the

## 798 SARS viral testing

tubing into the nostril. Aspirate nasopharyngeal secretions. Repeat this procedure for the other nostril. Collect the specimens in sterile vials.

- To obtain a *nasopharyngeal* or *oropharyngeal swab*, use only sterile Dacron or rayon swabs with plastic shafts. Do not use a cotton swab or swabs with wooden sticks because they may contain substances that inactivate some viruses and inhibit polymerase chain reaction testing. Insert the swab into the nostril. Leave the swab in place for a few seconds to absorb secretions. Swab both nostrils. (For *oropharyngeal culture*, swab the posterior pharynx and tonsillar areas, avoiding the tongue.)
- To collect *sputum*, educate the patient about the difference between sputum and oral secretions. Have the patient rinse the mouth with water and then expectorate deep-cough sputum directly into a sterile screw-cap sputum collection cup or sterile, dry container.
- To collect *blood*, collect whole blood in a serum separator tube. Collect blood in an EDTA (ethylenediaminetetraacetic acid; purple-top) tube for plasma testing.

### After

- Provide acute care for respiratory illness.
- If shipping the specimen domestically, use cold packs to keep the sample at 4° C. If shipping internationally, pack in dry ice.

## Abnormal findings

### SARS

### notes

**scrotal ultrasound** (Ultrasound of testes)

**Type of test** Ultrasound

**Normal findings**

Normal size, shape, and configuration of the testicles

**Test explanation and related physiology**

Present uses for scrotal ultrasound include:

- Evaluation of scrotal masses
- Measurement of testicular size
- Evaluation of scrotal trauma
- Evaluation of scrotal pain and identification of torsion of the testicle
- Evaluation of occult testicular neoplasm
- Surveillance of patients with prior primary or metastatic testicular neoplasms
- Follow-up for testicular infections
- Location of undescended testicles

The testicle and extratesticular intrascrotal tissues are examined. Both benign and malignant tumors can be identified with ultrasound. Benign abnormalities such as testicular abscess, orchitis, testicular infarction, and testicular torsion also can be identified. Extratesticular lesions such as hydrocele (fluid in the scrotum), hematocele (blood in the scrotum), and pyocele (pus in the scrotum) can be identified.

Ultrasonography of the scrotum is now the preferred method to identify torsion of the testicle. Ultrasonography is a very accurate method of identifying microlithiasis in the testicles. When identified, microcalcifications in the testicle indicate a marked increased risk for testicular cancer. Calcifications can also occur after orchitis or trauma. In most cases, both testicles are routinely imaged during the ultrasound examination.

The use of color Doppler is very helpful in determining blood flow to the testicle. If there is torsion of the testicle, color Doppler will indicate markedly reduced blood flow; immediate surgical exploration would be required. Scrotal ultrasound has replaced scrotal nuclear imaging for the diagnosis of testicular torsion because results can be obtained immediately.

There is very little discomfort associated with testicular ultrasound. It is usually performed by an ultrasound technologist and interpreted by an ultrasound physician.

## 800 scrotal ultrasound

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting is required.

#### During

- Note the following procedural steps:
  1. Careful examination of the scrotum is performed by the physician. Usually a short history is obtained.
  2. The scrotum is supported by a towel or cradled by the examiner's gloved hand.
  3. A greasy, conductive paste is applied to the scrotum. This paste enhances sound wave transmission and reception.
  4. Thorough scanning in the sagittal, transverse, and oblique projections is performed.
- The test takes approximately 20 to 30 minutes.

#### After

- Remove the coupling agent from the patient's scrotum.

### Abnormal findings

Benign testicular tumor  
Malignant testicular tumor  
Occult testicular tumor  
Testicular infection (orchitis)  
Hydrocele  
Hematocele  
Pyocele  
Varicocele  
Epididymitis  
Spermatocele  
Scrotal hernia  
Cryptorchidism  
Hematoma  
Testicular torsion

#### notes

## semen analysis (Sperm count, Sperm examination)

**Type of test** Fluid analysis

### Normal findings

Volume: 2-5 mL

Liquefaction time: 20-30 minutes after collection

Appearance: Normal

Motile/mL:  $\geq 10 \times 10^6$

Sperm/mL:  $\geq 20 \times 10^6$

Viscosity:  $\geq 3$

Agglutination:  $\geq 3$

Supravital:  $\geq 75\%$  live

Fructose: Positive

pH: 7.12-8

Sperm count (density):  $\geq 20$  million/mL

Sperm motility:  $\geq 50\%$  at 1 hour

Sperm morphology:  $> 30\%$  (Kruger criteria  $> 14\%$ ) normally shaped

### Test explanation and related physiology

Semen production depends on the function of the testicles; semen analysis is a measure of testicular function. Gonadotropin-releasing hormone (Gn-RH) stimulates the pituitary to produce follicle-stimulating hormone (FSH) and luteinizing hormone (LH, also called interstitial cell-stimulating hormone). The FSH stimulates the Sertoli cell growth to encourage sperm production. LH stimulates the Leydig cells to produce testosterone, which in turn stimulates the seminiferous tubules to produce sperm. Inadequate sperm production can be the result of primary gonadal failure (because of age, genetic cause [Klinefelter syndrome], infection, radiation, or surgical orchiectomy) or secondary gonadal failure (because of pituitary diseases). These forms of gonadal failure can be differentiated by measuring LH and FSH levels. In primary gonadal failure, LH and FSH levels are increased. In secondary gonadal failure, they are decreased. Stimulation tests using Gn-RH agonists such as leuprolide acetate clomiphene, or human chorionic gonadotropin are also used in the differentiation. Men with *aspermia* (no sperm) or *oligospermia* ( $< 20$  million/mL) should be evaluated endocrinologically for pituitary, thyroid, or testicular aberrations.

Semen analysis is one of the most important aspects of the fertility workup because the cause of a couple's inability to conceive often lies with the man. After 2 to 3 days of sexual abstinence,

semen is collected and examined for volume, sperm count, motility, and morphology.

The freshly collected semen is first measured for *volume, pH, and viscosity*. After liquefaction of the white, gelatinous ejaculate, a sperm count is done. Men with very low or very high counts likely are infertile. The motility of the sperm is then evaluated; at least 50% should show rapid ( $> 25 \mu\text{m/s}$  at  $37^\circ \text{C}$ ) or sluggish progressive motility. Morphology is studied by staining a semen preparation and calculating the number of sperm with normal versus abnormal morphology. Using the Kruger criteria, sperm morphology must be greater than 14% to be considered normal. Morphology of less than 4% is associated with severe infertility. More exhaustive semen analysis or second-tier testing for male infertility may include sperm functional testing, identification of sperm antibodies (see p. 93), and biochemical testing.

Sperm functional tests include:

1. Sperm–cervical mucus interaction: This is a postcoital test that evaluates the sperm–cervical mucus interaction. Normal is more than 10 to 20 motile sperm per high-power field (hpf).
2. Computer-assisted semen analysis: In this test, several different sperm kinetics are evaluated, including velocities, linearity, and amplitude of sperm head displacement.
3. Sperm penetration assay (SPA): This is a multistep laboratory test that offers a biologic assessment of several aspects of human sperm fertilizing ability.
4. Hemizona and zona pellucida binding tests: These include the hemizona assay (HZA) and a competitive intact zona binding assay. These tests evaluate the interaction between the spermatozoa and the zona pellucida of the female egg. These tests thereby evaluate many functions of the sperm at one time.

Interruption in sperm DNA integrity is a potential cause of male infertility. Although sperm with fragmented DNA may be able to fertilize oocytes, subsequent embryo and fetal development may be impaired. DNA fragmentation in sperm increases with age. Therefore, impaired DNA integrity may be an increasing infertility factor among older couples. Testing for DNA integrity include *sperm chromatin structure assay test* and the *sperm DNA fragmentation assay (SDEA) test*. The sperm specimen is considered abnormal if more than 70% of the sperm have abnormal forms. There are several other direct and indirect tests of DNA damage and chromosomal tests measuring chromosomal numerical abnormalities in sperm.

Sperm biochemical testing includes measurements of zinc, citric acid, glucosidase and the *Hyaluronan binding assay* (HBA). The HBA is based on the ability of mature, but not immature, sperm to bind to hyaluronan, the main mucopolysaccharide of the egg matrix and a component of human follicular fluid. Hyaluronan-binding capacity is acquired late in the sperm maturation process; immature sperm lack this ability. Therefore a low level of sperm binding to hyaluronan suggests that there is a low proportion of mature sperm in the sample. Similar to the sperm penetration assay, it has been suggested that the HBA assay may be used to determine the need for an intracytoplasmic sperm injection procedure as part of an assisted reproductive technique.

A single sperm analysis, especially if it indicates infertility, is inconclusive because sperm count varies from day to day. A semen analysis should be done at least twice and possibly a third time, 3 weeks apart. A normal semen analysis alone does not accurately assess the male factor unless the effect of the partner's cervical secretion on sperm survival is also determined. Sperm antibody testing (see p. 93) is also performed on the specimen.

In addition to its value in infertility workups, semen analysis is also helpful in documenting adequate sterilization after a vasectomy. It is usually performed 6 weeks after the surgery. If any sperm are seen, the adequacy of the vasectomy must be suspect.

### Interfering factors

- ✦ Drugs that may cause *decreased* semen levels include antineoplastic agents (e.g., methotrexate), cimetidine, estrogens, and methyltestosterone.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to abstain from sexual activity for 2 to 3 days before collecting the specimen. Prolonged abstinence before the collection should be discouraged because the quality of the sperm cells, especially their motility, may diminish.
  - Give the patient the proper container for the sperm collection.
- PT** Instruct the patient to avoid alcoholic beverages for several days before the collection.
  - For evaluation of the adequacy of vasectomy, the patient should ejaculate once or twice before the day of examination.

## 804 semen analysis

### During

- Note that semen is best collected by ejaculation into a clean container. For the best results, the specimen should be collected in the physician's office or laboratory by masturbation.
- PT Note that less satisfactory specimens can be obtained in the patient's home by coitus interruptus or masturbation. Note the following procedural steps:
  1. Instruct the patient to deliver these home specimens to the laboratory within 1 hour after collection.
  2. Tell the patient to avoid excessive heat and cold during transportation of the specimen.

### After

- Record the date of the previous semen emission along with the collection time and date of the fresh specimen.
- PT Tell the patient when and how to obtain the test results. Remember that abnormal results may have a devastating effect on the patient's sexuality.

### Abnormal findings

Infertility

Vasectomy

Orchitis

Testicular atrophy

Testicular failure

Hyperpyrexia

### notes



**sentinel lymph node biopsy (SLNB, Lymphoscintigraphy)**

**Type of test** Nuclear scan

**Normal findings**

Uptake is noted in one or more lymph nodes. No tumor in the sentinel node.

**Test explanation and related physiology**

Lymphoscintigraphy is used to identify the sentinel lymph node, which is the first lymph node in line to catch metastasis from a nearby primary tumor. It is primarily used in breast cancer and melanoma. With this procedure, the sentinel node is identified and biopsied. This test is an important part of the standard treatment for breast and melanoma cancer surgery.

**Contraindications**

- Patients who have a large cancer in which lymph node metastasis is very likely
- Patients in early pregnancy unless the benefit outweighs the risk of damage to the fetus

**Potential complications**

- Anaphylaxis with injection of isosulfan blue dye

**Procedure and patient care****Before**

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- Because this is an operative procedure, routine preoperative nursing processes should be carried out, including obtaining operative consent, keeping the patient NPO (nothing by mouth), and surgical site preparation as ordered.

**During**

- Note the following procedural steps:
  1. After the injection of technetium, the lymph node drainage basin is then scanned immediately and 1 to 24 hours later.
  2. In the operating room, a handheld gamma detector locates hot areas of radionuclide uptake in the lymph node-bearing area. The most proximal hot node is excised as the sentinel node.

## 806 sentinel lymph node biopsy

3. When using isosulfan blue, 4 to 5 mL of dye is injected around the tumor. Several minutes later, the lymphatics are stained blue and can be identified for biopsy.

### After

- PT** Inform the patient that no precautions are required if technetium is used because the radionuclide dose is minimal.
- PT** If isosulfan blue dye is used, the patient's skin may develop a transient blue hue (looking almost like severe cyanosis). This will dissipate over the next 6 hours.
- PT** Warn the patient that the urine will have a blue tinge as a result of the isosulfan blue dye injection.
  - Observe the patient for signs of allergy (rare) caused by the blue dye injection.

### Abnormal findings

Metastatic tumor to lymph node

### notes

## serotonin (5-hydroxytryptamine, 5-HT) and chromogranin A

**Type of test** Blood

### Normal findings

Chromogranin A:  $\leq 225$  ng/mL

Serotonin:  $\leq 230$  ng/mL

### Test explanation and related physiology

Serotonin is synthesized from the essential amino acid tryptophan chiefly in the gastrointestinal (GI) enterochromaffin (EC) cells. Many different stimuli can release serotonin from EC-cells. After it is secreted, in concert with other gut hormones, serotonin increases GI blood flow, motility, and fluid secretion. On the first pass through the liver, 30% to 80% of serotonin is metabolized, predominately to 5-hydroxyindoleacetic acid (5-HIAA) (see p. 528), which is then excreted by the kidneys.

The main diseases that may be associated with measurable increases in serotonin are neuroectodermal tumors, in particular tumors arising from EC cells. These tumors are collectively referred to as *carcinoids*. Most symptoms of carcinoid tumors are caused by elevated serotonins (carcinoid syndrome). The carcinoid syndrome consists of flushing, diarrhea, right-sided valvular heart lesions, and bronchoconstriction.

Diagnosis of carcinoid tumors with symptoms suggestive of carcinoid syndrome rests on measurements of serum serotonin, urinary 5-HIAA (see p. 528), and *serum chromogranin A* (a peptide that is cosecreted alongside serotonin by the neuroectodermal cells). Metastasizing midgut carcinoid tumors usually produce blood or serum serotonin concentrations greater than 1000 ng/mL. It is only after carcinoid tumors metastasize that serotonin becomes detectable.

Disease progression can be monitored in patients with serotonin-producing carcinoid tumors by measurement of serotonin or chromogranin A in blood.

Chromogranin A also acts as a useful cancer tumor marker (see p. 194) for other neuroendocrine neoplasms, including carcinoids, pheochromocytomas, neuroblastomas, medullary thyroid carcinomas, some pituitary tumors, functioning and nonfunctioning islet-cell tumors, and other amine precursor uptake and decarboxylation (APUD) tumors. It can also serve as a sensitive means for detecting residual or recurrent disease in treated patients. Carcinoid tumors, in particular colon and rectal carcinoids, almost always secrete chromogranin A.

### Interfering factors

- ▼ Drugs that may cause *increased* serotonin levels include lithium, monoamine oxidase inhibitors, methyldopa, morphine, and reserpine.
- ▼ Drugs that *decrease* serotonin levels include selective serotonin reuptake inhibitors (e.g., fluoxetine).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

- Carcinoid tumors
- Neuroendocrine tumors
- Pheochromocytoma
- Small cell lung cancer

### notes

## sexual assault testing (rape testing)

**Type of test** Blood; fluid analysis

### Normal findings

No physical evidence of sexual assault

### Test explanation and related physiology

The sexual assault victim needs to have psychoemotional support, treatment of any physical injuries, and accurate and reliable evidentiary testing. Nearly all acute care centers have protocols in place that provide that care to victims of sexual assault. Furthermore, in most circumstances, there are nurses specifically trained in obtaining the appropriate specimens. These nurses know the importance of following the chain of evidence protocols to ensure that evidence is admissible in court.

The patient is first interviewed in a nonjudgmental manner. A thorough gynecologic history is obtained. A brief summary of the assault (if there was vaginal, oral, or anal penetration) and timing of the assault is important. After 72 hours, very little evidence persists. It is important to ascertain if the victim changed clothing, showered, or used a douche before coming to the hospital. These will affect the presence of evidence. The general demeanor of the patient, status of the clothing, and physical maturation assessment is documented.

The victim's clothes are removed and separately placed in a paper bag for possible DNA sources of the victim's or assailant's body parts. Plastic bags are not used because bacteria may grow in them and can destroy DNA. Photographs of all injuries should be obtained, if possible. The victim is then examined for signs of external and internal injuries. A pelvic examination is then performed. A "sexual assault evidence collection kit" (also known as a rape kit, sexual assault kit [SAK], a sexual assault forensic evidence [SAFE] kit, a sexual assault evidence collection kit [SAECK], a sexual offense evidence collection [SOEC] kit, or a physical evidence recovery kit [PERK]) is used to obtain all the needed specimens. The directions must be carefully followed to ensure that any and all evidence is obtained and is useful toward identification and conviction of any perpetrator.

Vaginal secretions (or from any area of penetration) are obtained for sperm (see p. 801), or other cells from the assailant. Acid phosphatase (see p. 7) or prostate specific antigen (PSA) (see p. 743) are also obtained using this specimen. Cervical

## 810 sexual assault testing

secretions are obtained for sexually transmitted disease (STD) (see p. 813) testing. Wet preps may show motile sperm. These anatomic areas along with the anorectal area are swabbed per directions in the kit. In a male victim, penile and anorectal areas are swabbed. Pubic hair is obtained by combing or plucking. STD testing includes syphilis (see p. 859), trichomoniasis (see p. 813), gonorrhea (see p. 813), and chlamydia (see p. 230). Later, blood can be tested for human immunodeficiency virus (HIV) (see p. 506), hepatitis (see p. 493), and pregnancy (see p. 515 for women of childbearing age).

Next, blood specimens are obtained for DNA testing per the testing kit directions. More blood or urine may also be collected for evidence of mind-altering drugs (e.g., flunitrazepam, alcohol, or benzodiazepines). After this testing, a more detailed examination of the vagina, cervix, and rectum is performed using a Wood lamp to more easily identify saliva or sperm from the assailant. These areas are examined for subtle injuries from forced penetration. Two methods used to identify these injuries are the toluidine blue dye test and use of a colposcope (see p. 262). The *toluidine blue dye test* can also be used to identify recent or healed genital or anorectal injuries. A 1% aqueous solution is applied to the area of concern and washed off with a lubricant (e.g., K-Y Jelly) or a 1% acetic acid solution. Injured mucosa will retain the dye and become more apparent to the naked eye. Finally, the fingernails are scraped underneath because they may potentially contain tissue from the assailant. On completion of the examination, the victim is usually interviewed by the police for further investigation.

Unless medically contraindicated, all victims should be offered antimicrobial therapy to prevent STDs. The use of anti-retroviral drugs in the prevention of HIV transmission may be recommended, and the current guideline for postexposure prophylaxis after needle stick injuries should be used. It may also be advisable to offer victims a hepatitis B vaccination or hepatitis B immunoglobulin. Victims who are at risk for HIV infection and AIDS should be given counseling.

A pregnancy test should be done before any treatment or drugs are prescribed. If there is a risk of pregnancy, the victims should be offered postcoital contraception if the rape occurred less than 72 hours before examination by the health worker. If it occurred more than 72 hours but less than 7 days before the examination, an intrauterine contraceptive device may be used to prevent pregnancy. Pregnancy testing may be repeated in the succeeding week after the rape.

Clinicians may occasionally be called on to perform forensic evaluation of an assault perpetrator. Principles for evidentiary examinations are similar to those for victims and require evidence collection kits and strict attention to maintain the chain of evidence. Swabs, hair combing, and fingernail sampling are obtained. Penile swabs should be collected from the shaft, glans, and area under the foreskin; finger swabs should be done in cases of digital penetration of a victim. Bruises, scratches, and bite marks are identified, with swabbing of bites and scratches to identify victim DNA. Blood samples for HIV and hepatitis B can be drawn and held.

### Contraindications

- The patient is emotionally unable to undergo testing.

### Interfering factors

- Delays in examination after the alleged attack diminish the possibility of identifying meaningful evidence.

### Procedure and patient care

#### Before

**PT** Explain the procedure and provide emotional support.

- Obtain consent to treat the patient.
- Assess the patient's emotional condition and determine whether the victim is able to undergo sexual assault testing.

#### During

- Use the SAECK or similar test kit exactly as described to maintain the chain of evidence (see [Box 6](#)).
- Properly handle specimens to maintain the chain of custody.

#### **BOX 6 DNA evidence collection: special precautions**

To avoid contamination of evidence that may contain DNA, the special sexual assault kit should be used and the following precautions taken:

- Wear gloves and change them often.
- Use disposable instruments or clean them thoroughly before and after handling each sample.
- Avoid touching any area where DNA may be present.
- Avoid talking, sneezing, or coughing over evidence.
- Avoid touching your face, nose, and mouth when collecting and packaging evidence.
- Keep evidence such as clothing and underwear dry and transport it at room temperature.
- Ensure that the chain of custody is maintained at all times.

## 812 sexual assault testing

- Refrigerate all blood and urine samples containing biological evidentiary material such as DNA to prevent putrefaction.
- It is important to examine carefully all areas of the body to help corroborate the victim's version of the alleged events.

### After

- Notify police of the alleged assault.
- Assess the patient's need for urgent counseling support and make arrangements as needed.

### Abnormal findings

Rape

Sexual assault

notes



## sexually transmitted disease testing (STD testing)

**Type of test** Microscopic examination

### Normal findings

No evidence of STD

### Test explanation and related physiology

In the United States, common STDs include *Chlamydia*, genital herpes (herpes simplex virus), human papilloma virus (HPV), syphilis, human immunodeficiency virus (HIV), trichomonas, and gonorrhea (Table 35). In this test discussion, we will concentrate on *Trichomonas vaginalis* and *Neisseria gonorrhoeae*, as all others are discussed elsewhere in this reference book. Early identification of STDs enables sexual partners to obtain treatment as soon as possible and thereby reduce the risk of disease spread.

**TABLE 35** Sexually transmitted diseases (STDs) and methods of diagnosis

Disease	Method of diagnosis
Gonorrhea	Cervical, urethral, anal, oropharyngeal cultures
<i>Chlamydia</i>	Cervical and urethral cultures, serology,
<i>Lymphogranuloma venereum</i>	DNA probe testing (see p. 230)
<i>C. trachomatis</i>	
Herpes genitalis	Culture from lesion, serology (see p. 498)
Syphilis	Serology, fluid cultures (central nervous system), darkfield slide (see p. 859)
Hepatitis	Serology, nucleic acid testing (see p. 493)
HIV	Serologic, virologic, nucleic acid testing (see p. 503, 506)
<i>Trichomonas vaginalis</i>	Cervical and urethral cultures, urine, ThinPrep Pap, serology, nucleic acid amplification tests
<i>Candida</i>	Wet mount, fungal culture
<i>Gardnerella vaginalis</i>	Cervical, urethral, anal cultures

## 814 sexually transmitted disease testing

Furthermore, prompt treatment reduces the risk of infertility in women. If the STD result is positive, sexual partners should be evaluated and treated. Performing STD testing is also part of the prenatal workup.

*T. vaginalis* can cause urethritis, vaginitis, endometritis, pelvic inflammatory disease, pharyngitis, proctitis, epididymitis, prostatitis, and salpingitis. Children born of infected mothers may develop conjunctivitis, pneumonia, neonatal blindness, or neonatal neurologic injury and may even die. The most commonly used method for detection is *microscopic examination of a wet-mount preparation of vaginal secretions*. However, this method has only a 35% to 80% sensitivity. *Culture of urethral or vaginal secretions* also suffers from relatively low sensitivity. Culture is technically challenging and takes 5 to 7 days to complete. Molecular methods of testing urethral and vaginal secretions offer high sensitivity and specificity for detection of trichomoniasis.

Gonorrhea is caused by the bacterium *N. gonorrhoeae*. Many infections in women are asymptomatic. This organism causes genitourinary infections in women. Because infection in men is commonly associated with symptoms, screening of asymptomatic patients is not indicated. However, in light of the risk for asymptomatic infection in women, screening is recommended for women at high risk for infection. High-risk women include women with previous gonorrhea or other STD, inconsistent condom use, and new or multiple sex partners and women in certain demographic groups such as those in communities with high STD prevalence.

*Culture* was previously considered to be the gold standard test for diagnosis of *N. gonorrhoeae* infection. Yet successful culture methods are difficult. Molecular laboratory methods performed on urethral, rectal, vaginal, or oropharyngeal secretions provide superior sensitivity and specificity.

To obtain an appropriate specimen for women, swabs (that are sometimes specific to the particular laboratory) are obtained from the endocervix, vagina, urethra, urine, or a Pap ThinPrep. For men, a swab of the urethra or a urine specimen is used for testing. Rectal and throat swabs are performed in persons who have engaged in anal and oral intercourse. Because rectal gonorrhea accompanies genital gonorrhea in a high percentage of women, rectal cultures are recommended in all women with suspected gonorrhea. If the STD culture result is positive, treatment during pregnancy can prevent possible fetal complications (e.g., ophthalmia neonatorum) and maternal complications. Rectal and orogastric specimens should be performed on the neonates of infected mothers.

STD cultures and smears are obtained by a physician or nurse in several minutes during a pelvic examination. Very little discomfort is associated with these procedures.

### Interfering factors

- *N. gonorrhoeae* is very sensitive to lubricants and disinfectants.
- Menses may alter test results.
- In women, douching within 24 hours before a cervical culture makes fewer organisms available for culture.
- In men, voiding within 1 hour before a urethral culture washes secretions out of the urethra.
- Fecal material may contaminate an anal culture.

### Procedure and patient care

#### Before

- PT** Explain the purpose and procedure to the patient. Use a matter-of-fact, nonjudgmental approach.
- PT** Tell the patient that no fasting or sedation is required.

#### During

##### *Cervical culture*


- The female patient is told to refrain from douching and tub bathing before the cervical culture.
- The patient is placed in the lithotomy position, and a moistened, unlubricated vaginal speculum is inserted to expose the cervix (see [Figure 31](#), p. 669).
- Excess cervical mucus is removed with a cotton ball.
- A sterile cotton-tipped swab is inserted into the endocervical canal and moved from side to side to obtain the specimen.
- The swab is placed in sterile saline or a transporting fluid obtained from the laboratory. The specimen should be plated as soon as possible.

##### *Anal canal culture*

- An anal culture of the female or male patient is taken by inserting a sterile, cotton-tipped swab approximately 1 inch into the anal canal ([Figure 38](#)).
- If stool contaminates the swab, a repeat swab is taken.

##### *Oropharyngeal culture*

- This culture should be obtained in male and female patients who have engaged in oral intercourse.
- A throat culture is best obtained by depressing the patient's tongue with a wooden tongue blade and touching the posterior wall of the throat with a sterile cotton-tipped swab.

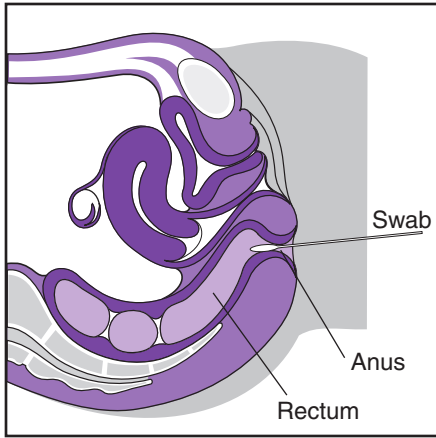

**FIG. 38** Rectal culture of the female. Method for obtaining an anorectal culture for sexually transmitted diseases in a female patient.

*Urethral culture*


- The urethral specimen should be obtained from the male patient before he voids. Voiding within 1 hour before collection washes secretions out of the urethra, making fewer organisms available for culture. The best time to obtain the specimen is before the first morning micturition.
- A culture is taken by inserting a sterile swab gently into the anterior urethra (Figure 39).
- It is advisable to place the male patient in the supine position to prevent falling if vasovagal syncope occurs during introduction of the cotton swab or wire loop into the urethra.
- The patient is observed for hypotension, bradycardia, pallor, sweating, nausea, and weakness.
- In a male patient, prostatic massage may increase the chances of obtaining positive cultures.

*Urine culture*

- Obtain the first catch voided specimen in the female patient. (Urine cultures for STD are not helpful in male patients.)

*Pap smear ThinPrep*

See p. 667.

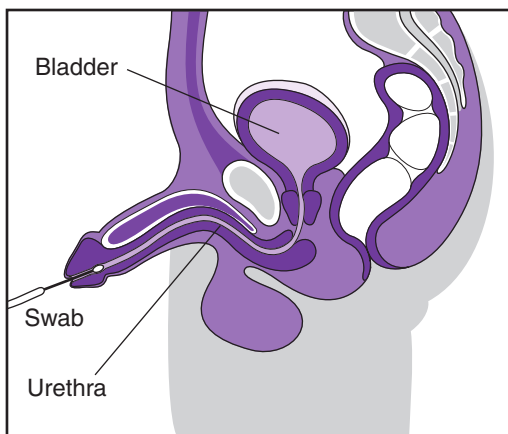

**FIG. 39** Urethral culture of the male. Method for obtaining a urethral culture for sexually transmitted diseases in a male patient.

### After

- Place the swabs for gonorrhea in a Thayer-Martin medium and roll them from side to side.
  - Label and send the culture bottle to the microbiology laboratory.
  - Transport the specimen to the laboratory as soon as possible.
  - Handle all specimens as though they were capable of transmitting disease.
  - Do not refrigerate the specimen.
  - Mark the laboratory slip with the collection time, date, source of specimen, patient's age, current antibiotic therapy, and clinical diagnosis.
- PT** Advise the patient to avoid intercourse and all sexual contact until test results are available.
- PT** If the culture results are positive, tell the patient to receive treatment and to have sexual partners evaluated.
- Note that repeat cultures should be taken after completion of treatment to evaluate therapy.

### Abnormal findings

STDs

notes

## sialography

**Type of test** X-ray

### Normal findings

No evidence of pathology in the salivary ducts and related structures

### Test explanation and related physiology

Sialography is an x-ray procedure used to examine the salivary ducts (parotid, submaxillary, submandibular, and sublingual) and related glandular structures after injection of a contrast medium into the desired duct. This procedure is used to detect calculi, strictures, tumors, or inflammatory disease in patients who complain of pain, tenderness, or swelling in these areas.

### Contraindications

- Patients with mouth infections

### Potential complications

- Allergic reaction to the iodinated dye. This rarely occurs because the dye is not administered intravenously.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. The thought of a dye injection in the mouth is frightening to many patients. Provide emotional support. See p. xviii for radiation exposure and risks.
- Obtain informed consent if required by the institution.
- PT** Instruct the patient to remove jewelry, hairpins, and dentures, which could obscure x-ray visualization.
- PT** Instruct the patient to rinse his or her mouth with an antiseptic solution before the procedure to reduce the possibility of introducing bacteria into the ductal structures.

#### During

- Note the following procedural steps:
  1. X-ray studies are taken before the dye injection to ensure that stones are not present, which could prevent the contrast material from entering the ducts.
  2. The patient is placed in a supine position on an x-ray table.
  3. The contrast medium is injected directly into the desired orifice via a cannula or a special catheter.
  4. X-ray images are taken with the patient in various positions.

5. The patient is given a sour substance (e.g., lemon juice) orally to stimulate salivary excretion.
  6. Another set of x-ray studies is taken to evaluate ductal drainage.
- Note that a radiologist performs this procedure in the radiology department in less than 30 minutes.
- PT** Tell the patient that he or she may feel a little pressure as the contrast medium is injected into the ducts.

### After

- PT** Encourage the patient to drink fluids to eliminate the dye.

### Abnormal findings

Calculi  
Strictures  
Tumor  
Inflammatory disease

### notes

**sickle cell screen** (Sickledex, Hemoglobin [Hgb] S test)**Type of test** Blood**Normal findings**

No sickle cells present or no Hgb S identified

**Test explanation and related physiology**

Both sickle cell disease (homozygous for Hgb S) and sickle cell trait (heterozygous for Hgb S) can be detected by this screening study. Sickle cell anemia results from a genetic homozygous defect and is caused by the presence of Hgb S instead of Hgb A. When Hgb S becomes deoxygenated, it tends to bend in a way that causes the red blood cell (RBC) to assume a sickle shape. Hgb S is found in varying quantities in 8% to 10% of the black population.

The *Sickledex test* is only a screening test, and its sensitivity varies according to the method used by the laboratory. Double heterozygosity for sickle trait when combined with another hemoglobinopathy (e.g., Hgb C disease) can cause a sickling disease. The definitive diagnosis of sickle cell disease or trait is made by Hgb electrophoresis (see p. 490) in which Hgb S can be identified and quantified.

**Interfering factors**

- Any blood transfusions within 3 to 4 months before the sickle cell test may cause false-negative results because the donor's normal hemoglobin may dilute the recipient's Hgb S.
- Polycythemia or paraproteinemias may cause false-positive solubility results.
- Infants younger than 3 months may have false-negative results.
- ✚ Phenothiazines may cause false-negative results.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender
- If the test is positive, further testing is done.
- PT Inform patients with sickle cell anemia that they should avoid situations in which hypoxia may occur (e.g., strenuous exercise, travel to high-altitude regions).



## Abnormal findings

Sickle cell trait

Sickle cell anemia

## notes

## Sims-Huhner test (Postcoital test, Postcoital cervical mucus test, Cervical mucus sperm penetration test)

**Type of test** Fluid analysis

### Normal findings

Cervical mucus adequate for sperm transmission, survival, and penetration

6-20 active sperm per hpf

### Test explanation and related physiology

The Sims-Huhner test consists of a postcoital examination of the cervical mucus to measure the ability of the sperm to penetrate the mucus and maintain motility. This test is used in the evaluation of infertile couples. It evaluates interaction between the sperm and the cervical mucus. It also measures the quality of the cervical mucus. This test can determine the effect of vaginal and cervical secretions on the activity of the sperm. It is performed only after a previously performed semen analysis has been determined to be normal.

This test is performed during the middle of the ovulatory cycle because at this time, the secretions should be optimal for sperm penetration and survival. During ovulation, the quantity of cervical mucus is maximal, but the viscosity is minimal, thus facilitating sperm penetration. The *endocervical mucus sample* is examined for color, viscosity, and tenacity (*spinnbarkeit*). The fresh specimen is then spread on a clean glass slide and examined for the presence of sperm. Estimates of the total number of sperm and of the number of motile sperm per high-power field are reported. Normally 6 to 20 active sperm cells should be seen in each microscopic high-power field; if sperm are present but not active, the cervical environment is unsuitable (e.g., abnormal pH) for their survival.

After the specimen has dried on the glass slide, the mucus can be examined for *ferning*. This pattern is correlated with estrogen activity and is therefore present in all ovulatory women at midcycle. When the cervical mucus is checked again immediately before menstruation, no ferning is found because of progesterone activity.

This test has limited diagnostic potential and poor predictive value. Its use has been associated with increased testing without improvement in pregnancy rates. Furthermore, cervical factor infertility is easily addressed by performing intrauterine inseminations.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Inform the patient that basal body temperature recordings should be used to indicate ovulation.
- PT** Tell the patient that no vaginal lubrication, douching, or bathing is permitted until after the vaginal cervical examination because these factors will alter the cervical mucus.
- PT** Inform the patient that this study should be performed after 3 days of sexual abstinence.
- PT** Instruct the patient to remain in bed for 10 to 15 minutes after coitus to ensure cervical exposure to the semen. After resting, the patient should report to her physician for examination of her cervical mucus within 2 hours after coitus.

### During

- Note that with the patient in the lithotomy position, the cervix is exposed by an unlubricated speculum. The specimen is aspirated from the endocervix and delivered to the laboratory for analysis.
- Note that this procedure is performed by a physician in approximately 5 minutes.
- PT** Tell the patient that the only discomfort associated with this study is insertion of the speculum.

### After

- PT** Tell the patient how and when she may obtain the test results.

## Abnormal findings

Infertility

Suspected rape

notes

**skin biopsy** (Cutaneous immunofluorescence biopsy, Skin biopsy antibodies, Skin immunohistopathology, Direct immunofluorescence antibody test)

**Type of test** Microscopic examination

### Normal findings

Normal skin histology

No evidence of immunoglobulin (Ig) G, IgA, or IgM antibody or of complement C3 or fibrinogen

### Test explanation and related physiology

Autoimmune skin diseases are associated with autoantibodies in the skin and serum. Direct (testing for antibodies in the skin) immunofluorescence antibody (IFA) is most specific and diagnostic. For this study, a tissue specimen in or around the skin or mucosal lesion is obtained and evaluated by routine histology and by IFA methods for deposition of human immunoglobulins (IgG, IgA, or IgM), complement C3, or fibrinogen components.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Obtain an informed consent.

#### During

- The skin area used for biopsy is surgically biopsied.
- A 4-mm punch biopsy or elliptical tissue excision is obtained.

#### After

- Apply a dry, sterile dressing over the biopsy site.

**PT** Tell the patient that results may not be available for several days.

- Deliver the specimen to the laboratory immediately after the biopsy is taken.

### Abnormal findings

Systemic lupus erythematosus

Discoid lupus erythematosus

Pemphigus

Bullous pemphigoid

Dermatitis herpetiformis

#### notes

## skull x-ray

**Type of test** X-ray

### Normal findings

Normal skull and surrounding structures

### Test explanation and related physiology

An x-ray image of the skull allows for visualization of the bones making up the skull, the nasal sinuses, and any cerebral calcification. Skull x-rays are rarely indicated today because of the availability of CT scanning of the brain (see p. 278). However, skull x-rays are still used for determining skull bone suture lines in the evaluation of children with abnormal head shape or size.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Instruct the patient to remove all objects above the neck because metal objects and dentures prevent x-ray visualization of the structures they cover.
  - Avoid hyperextension and manipulation of the head if surgical injuries are suspected.
- PT** Tell the patient that no sedation or fasting is required.

#### During

- Note that the patient is taken to the radiology department and placed on an x-ray table. Axial, half-axial, posteroanterior, and lateral views of the skull are usually taken.
- Note that a technologist takes skull x-rays in a few minutes.
- Tell the patient that this test is painless.

#### After

- If a prosthetic eye is present, note this on the x-ray examination request.

### Abnormal findings

Skull fracture

Metastatic tumor

Sinusitis

Hemorrhage

Tumor

Hematoma

Congenital anomaly

Paget disease of bone

#### notes

## sleep studies (Polysomnography [PSG])

**Type of test** Electrodiagnostic; various

### Normal findings

Respiratory disturbance index (RDI): < 5 episodes of apnea per hour

Normal progress through sleep stages

No interruption in nasal or oral airflow

End-tidal CO<sub>2</sub>: 30-45 mm Hg

Oximetry: ≥ 90%; no oxygen desaturation of > 5%

Minimal snoring sounds

Electrocardiogram (ECG): no disturbances in rate or rhythm

No evidence of restlessness

No apnea

Multiple sleep latency test (MSLT): onset of sleep > 9 minutes

### Test explanation and related physiology

There are many types of sleep disorders. Sleep studies can identify the cause of the sleep disorders and indicate appropriate treatment. Sleep studies include polysomnography (PSG) and testing for wakefulness and sleepiness. A full PSG would include:

- *Electroencephalography*: This is limited to two or more channels (see p. 347).
- *Electrooculography*: This documents eye movements (see electroystagmography, p. 357).
- *Electromyography (EMG)*: This demonstrates muscle movement, usually of the chin and legs (see p. 350).
- *Electrocardiography*: This determines heart rate and rhythm (see p. 342).
- *Chest impedance*: This monitors chest wall movement and respirations.
- *Airflow monitors*: This measures the amount of airflow in and out of the mouth and nose.
- *CO<sub>2</sub> monitor*: This measures expiratory CO<sub>2</sub> levels.
- *Pulse oximetry*: This monitors tissue oxygen levels (see p. 658).
- *Sound sensors*: These are used to document snoring sounds.
- *Audiovisual recordings*: These are used to document restless motion and fitfulness.
- *Esophageal pH probe*: This is used only if gastroesophageal reflux is considered to be a cause of paroxysmal nocturnal dyspnea and coughing (see p. 377).

On occasions when sleep apnea alone is suspected, a four-channel PSG is performed. This more simplified test includes

the ECG, chest impedance, airflow monitor, and O<sub>2</sub> oximetry. Audiovisual recordings are also performed.

A sleep screening study is often performed to see whether full sleep studies are indicated. This is done during sleep using pulse oximetry. If no hypoxia occurs, significant sleep apnea would be rare, and full studies are not indicated.

Sleep apnea can be obstructive or central. Obstructive apnea is by far the most common and is caused by muscle relaxation of the posterior pharyngeal muscles. Breathing stops for 10 to 40 seconds. Central sleep apnea is highlighted by simple cessation of breathing not caused by an obstructed airway. Primary cardiac events that lead to significant and transient reduction in cardiac output can also cause apnea. Apnea from either cause is associated with an increase in heart rate, decreased oxygen levels, change in brain waves, and increased expiratory CO<sub>2</sub>. Obstructive apnea is also associated with progressively diminished airflow.

During a sleep study, electrodes for the ECG, electroencephalography (EEG), electrooculography, and EMG are applied. The chest impedance belt monitors are also placed. Under audiovisual monitoring, the patient is placed in a comfortable room and sleeps. During sleep, information is synchronously gathered. The EEG determines the various stages of sleep, and the physiologic changes during each stage are documented.

The study is usually completed in one night, although occasionally two nights are required. A second day is often required to administer the *multiple sleep latency test (MSLT)* or the *multiple wake test (MWT)*. The MSLT is a measure of the patient's ability to sleep during a series of structured naps. The MSLT is typically done in the morning. The MWT is a measure of the patient's ability to not fall asleep during a period of what should be wakefulness. These tests are used to diagnose narcolepsy that follows a night of inadequate sleep.

These tests can also be used to determine the success of therapy for sleep disorders. The sleep study may be repeated after the patient has started using continuous positive airway pressure (CPAP) or a dental fixture for therapy. When on therapy, no sleep apnea should be noted. If the sleep apnea is significant on the initial sleep study, a *split study* can be performed, in which sleep is evaluated while the patient uses a CPAP machine for the next 4 hours. During that time, appropriate CPAP settings are calibrated to reduce apneic episodes and minimize uncomfortable side effects.

Because of the expense and the psychoemotional difficulties associated with testing in a sleep laboratory, there has been

significant growth in unattended *home sleep studies*. The patient is attached to a multichannel monitor by a sleep technician as previously described. The technician does not remain in attendance. The monitoring device records all key data so that a sleep disorder can be identified.

*Actigraphy* can be used to determine sleep patterns and circadian rhythms. A sleep actigraph is a simple device that is worn like a wrist watch. It can be used during normal activities (except swimming or bathing) for several days and nights. It does not require an overnight stay at a sleep left. Doctors can use actigraphy to help diagnose sleep disorders, including circadian rhythm disorders such as jet lag and shift work disorders. This test can also detect how well sleep treatments are working. Actigraphy can be used with PSG or alone.

### **Interfering factors**

- Psychologically induced insomnia associated with being in a sleep left
- Environmental noises, temperature changes, or other sensations may affect the sleep pattern.
- Times for sleep testing that are different from usual times may affect sleep patterns and should be avoided.

### **Procedure and patient care**

#### **Before**

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to avoid caffeine for several days before testing.
- PT** Reassure the patient that monitoring equipment will not interrupt the patient's sleep pattern.
- Allow the patient to express concerns about videotaping and other forms of monitoring.
- Several sleep rating questionnaires are completed by both the patient and his or her sleeping partner.
- Age, weight, and medical history are recorded.

#### **During**

- Electrodes for ECG, EEG, and EMG are applied to the patient. Excess hair may need to be shaved on male patients.
- Airflow, oximetry, and impedance monitors are also applied.
- When the patient is comfortable, he or she is allowed to sleep.
- The lights are turned off, and monitoring begins.
- For PSG, the patient is asked to sleep per normal routine.



***Multiple sleep latency testing***

- The patient is asked to nap about every 2 hours throughout the testing period.
- The nap is terminated after 20 minutes.
- Between naps, the patient must stay awake.

***Multiple wake testing***

- The patient is asked to stay awake and not nap.
- Monitoring is similar to that described for PSG except for impedance, sound, and airflow monitors.

**After**

- On completion of the sleep cycle, the monitors and electrodes are removed.
- Test results take several days to collate and interpret.

**Abnormal findings**

Obstructive sleep apnea

Central sleep apnea

Cardiac sleep apnea

Insomnia

Narcolepsy

Restless legs syndrome

Parasomnia

REM disorder

**notes**

**small bowel follow-through** (SBF, Small bowel enema)

**Type of test** X-ray with contrast dye

**Normal findings**

Normal positioning, motility, and patency of the small intestine  
No evidence of intrinsic obstruction or extrinsic compression

**Test explanation and related physiology**

The SBF study is performed to identify abnormalities in the small bowel. Usually the patient is asked to drink barium; in patients who cannot drink, barium can be injected through a nasogastric tube. X-ray images are then taken at timed intervals (usually 30 minutes) to follow the progression of barium through the small intestine. Significant delays in transit time of the barium may occur with both benign and malignant forms of obstruction or diminished intestinal motility (ileus). On the other hand, the flow of barium is faster in patients who have hypermotility states of the small bowel (malabsorption syndromes). Failure of the progression through the small bowel can be seen in patients with partial mechanical small bowel obstruction or diminished intestinal motility, as seen in patients with diabetes. Furthermore, SBF series are helpful in identifying and defining the anatomy of small bowel fistulas.

A more accurate radiographic evaluation of the small intestine is provided by the *small bowel enema*. Unlike the SBF, in which the barium is swallowed by the patient, during the small bowel enema, the barium is injected into a tube previously passed to the small bowel. This small bowel enema provides better visualization of the entire small bowel because the barium is not diluted by gastric and duodenal juices. Tumors, ulcers, and small bowel fistulas are more easily identified and defined with the enema.

**Contraindications**

- Patients with a complete small bowel obstruction
- Patients suspected of having a perforated viscus
- Barium should not be used in these patients because it may cause prolonged and recurrent abscesses if it leaks out of the bowel. Gastrografin, a water-soluble contrast medium, can be used if perforation is suspected.
- Patients with unstable vital signs

**Potential complications**

- Barium-induced small bowel obstruction

## Interfering factors

- Barium in the intestinal tract from a previous barium x-ray image may obstruct adequate visualization of the small bowel.
- Food or fluid in the GI tract

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Instruct the patient not to eat anything for at least 8 hours before the test.
- PT** Inform the patient that the SBF series may take several hours. Suggest that the patient bring reading material or some paperwork to occupy his or her time.

### During

- Note the following procedural steps:
  1. A specially prepared drink containing barium sulfate is mixed as a milkshake, which the patient drinks through a straw.
  2. Usually, an upper GI series is performed concomitantly (see p. 922).
  3. The barium flow is followed through the upper GI tract fluoroscopically.
  4. At frequent intervals (15-60 minutes), repeat x-ray images are taken to follow the flow of barium through the small intestine. These images are repeated until barium is seen flowing into the right colon. This usually takes 60 to 120 minutes, but in patients with delayed progression of the barium, the test may take as long as 24 hours to complete.

### *Small bowel enema*

1. This is usually performed by placing a long weighted tube transorally; however, a tube also can be placed into the upper small bowel endoscopically.
  2. After the tube is in place, a thickened barium mixture is injected through the tube, and x-ray images are serially performed as described for the SBF.
- Note that this procedure is performed by a radiologist in the radiology department in approximately 30 minutes.
- PT** Tell the patient that this test is not uncomfortable.

### After

- PT** Inform the patient of the need to evacuate adequately all the barium. Cathartics (e.g., magnesium citrate) are recommended. Initially, stools will be white and should return to normal color with complete evacuation.

## 832 small bowel follow-through

### Abnormal findings

Small bowel tumor

Small bowel obstruction from intrinsic tumors

Small bowel obstruction from adhesions, extrinsic tumors, or hernia

Inflammatory small bowel disease (e.g., Crohn disease)

Malabsorption syndromes (e.g., Whipple disease, sprue)

Congenital anatomic anomaly (e.g., malrotation)

Congenital abnormalities (e.g., small bowel atresia, duplication, Meckel diverticulum)

Small bowel intussusception

Small bowel perforation

### notes

## small intestinal bacterial overgrowth test (SIBO test, Lactulose Breath Test)

**Type of test** miscellaneous

### Normal findings

Change in hydrogen < 20 ppm

Methane peak > 10 ppm

### Test explanation and related physiology

Small intestinal bacterial overgrowth (SIBO) is characterized by excessive bacteria or the presence of atypical microbiota in the small intestine. This breath test relies on measurement of gases produced by bacteria in the intestine—hydrogen ( $H_2$ ) and methane ( $CH_4$ )—after ingestion of lactulose in a fasting state. Lactulose is not absorbed by the intestines. Hydrogen gas is exclusively produced by colonic bacteria as a result of fermentation of the lactulose. In SIBO, the abnormal presence of abundant bacteria in the small intestine ferment the malabsorbed lactulose resulting in elevated concentration of exhaled hydrogen. An abnormal rise of hydrogen in the first 90 minutes of testing is positive for SIBO. Elevated exhaled methane is abnormal and is common in patients with irritable bowel disease, chronic constipation, or overeating.  $CO_2$  is routinely measured to determine sample integrity.

Conditions commonly associated with SIBO include irritable bowel syndromes, inflammatory bowel disease, celiac disease, diabetes, and obesity. These patients commonly complain of bloating, gas, pain, diarrhea, constipation, muscle pain (fibromyalgia), and bad breath. Bacterial overgrowth can also be identified by analyzing small intestinal fluid samples during enteroscopy.

In preparation for this test, the patient must wait at least 2 to 4 weeks after colonoscopy or barium enema. No antibiotics, bismuth containing agents, probiotics, stool softeners, laxatives, or antacids should be administered in that same time frame.

### Interfering factors

- Patients with lactose intolerance or diseases of maldigestion may have increased levels.
- ✚ The use of antibiotics, laxatives, antacids, or probiotics may decrease levels.

### Before

**PT** Explain the procedure to the patient. Inform the patient that multiple breath samples will be needed.

### **834 small intestinal bacterial overgrowth test**

**PT** Instruct the patient to fast for 12 hours before testing.

**PT** Tell the patient to avoid smoking, including second-hand smoke, for at least 1 hour before or at any time during the breath test.

#### **During**

- Provide a specified dose of lactulose usually diluted in 8oz of water.
- Collect breath samples at 20 minute intervals up to 60 to 120 minutes.

#### **After**

- Indicate actual time of the breath sample and the amount of time after ingestion of the sample.
- Immediately send test specimens to the laboratory.

#### **Abnormal findings**

Small bowel bacterial overgrowth

## sodium (Na), blood

**Type of test** Blood

### Normal findings

Adult/elderly: 136-145 mEq/L or 136-145 mmol/L (SI units)

Child: 136-145 mEq/L

Infant: 134-150 mEq/L

Newborn: 134-144 mEq/L

### Possible critical values

< 120 or > 160 mEq/L

### Test explanation and related physiology

Sodium is the major cation in the extracellular space, in which serum levels of approximately 140 mEq/L exist. Therefore sodium salts are the major determinants of extracellular osmolality. The sodium content of the blood is a result of a balance between dietary sodium intake and renal excretion.

Many factors regulate homeostatic sodium balance. Aldosterone causes conservation of sodium by decreasing renal losses. Natriuretic hormone, or third factor, increases renal losses of sodium. Antidiuretic hormone (ADH), which controls the resorption of water at the distal tubules of the kidney, also affects serum sodium levels.

Physiologically, water and sodium are very closely interrelated. As free body water is increased, serum sodium is diluted, and the concentration may decrease. The kidney compensates by conserving sodium and excreting water. If free body water were to decrease, the serum sodium concentration would rise; the kidney would then respond by conserving free water.

### Interfering factors

- Recent trauma, surgery, or shock may cause increased levels.
- ✚ Drugs that may cause *increased* levels include anabolic steroids, antibiotics, carbenicillin, clonidine, corticosteroids, cough medicines, estrogens, laxatives, methyldopa, and oral contraceptives.
- ✚ Drugs that may cause *decreased* levels include angiotensin-converting enzyme inhibitors, captopril, carbamazepine, diuretics, haloperidol, heparin, nonsteroidal antiinflammatory drugs, intravenous (IV) fluids, sulfonyleureas, triamterene, tricyclic antidepressants, and vasopressin.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or green

## Abnormal findings

### ▲ Increased levels (hypernatremia)

#### *Increased sodium intake*

Excessive dietary intake  
Excessive sodium in IV fluids

#### *Decreased sodium loss*

Cushing syndrome  
Hyperaldosteronism

#### *Excessive free body water loss*

Excessive sweating  
Extensive thermal burns  
Diabetes insipidus  
Osmotic diuresis  
GI loss

### ▼ Decreased levels (hyponatremia)

#### *Decreased sodium intake*

Deficient dietary intake  
Deficient sodium in IV fluids

#### *Increased sodium loss*

Addison disease  
Diarrhea  
Vomiting or nasogastric aspiration  
Diuretic administration  
Chronic renal insufficiency

#### *Increased free body water*

Excessive oral water intake  
Excessive IV water intake  
Congestive heart failure  
Syndrome of inappropriate ADH (SIADH) secretion  
Osmotic dilution

#### *Third-space losses of sodium*

Ascites  
Peripheral edema  
Pleural effusion  
Intraluminal bowel loss (ileus or mechanical obstruction)



## sodium (Na), urine

**Type of test** Urine (24-hour)

### Normal findings

40-220 mEq/day or 40-220 mmol/day (SI units)

Spot urine: > 20 mEq/L

Fractional excretion of sodium ( $FE_{Na}$ ): 1%-2%

### Test explanation and related physiology

This test evaluates sodium balance in the body by determining the amount of sodium excreted in urine over 24 hours. Sodium is the major cation in the extracellular space. Measuring the amount of sodium in the urine is useful for evaluating patients with volume depletion, acute renal failure, adrenal disturbances, and acid-base imbalances. In the setting of acute renal failure, whereas an increased value will indicate acute tubular necrosis, a low value would be typical of prerenal azotemia.

This test is also useful when the serum sodium concentration is low. For example, in patients with hyponatremia caused by inadequate sodium intake, urine sodium will be low. In patients with hyponatremia caused by chronic renal failure, however, urine sodium concentration will be high.

Urine sodium excretions are helpful when the urine output is low (< 500 mL/24 hr). However, a more accurate test to determine the cause of reduced urine output is the  $FE_{Na}$ .  $FE_{Na}$  is the fraction of sodium actually excreted relative to the amount filtered by the kidney.  $FE_{Na}$  is a calculation based on the concentrations of sodium (Na) and creatinine (Cr) in the blood and the urine as follows:

$$\text{Fractional excretion of sodium } (FE_{Na}) = (U_{Na} \times P_{Cr}) / (P_{Na} \times U_{Cr}) \times 100$$

$FE_{Na}$  is usually greater than 3% with acute tubular necrosis and severe obstruction of the urinary drainage of both kidneys. It is generally less than 1% in patients with acute glomerulonephritis, hepatorenal syndrome, and states of prerenal azotemia (e.g., congestive heart failure and dehydration).  $FE_{Na}$  may also be less than 1% with acute partial urinary tract obstruction.

### Interfering factors

- Dietary salt intake may increase sodium levels.
- Altered kidney function may affect levels.

## 838 sodium

- ✚ Drugs that may cause *increased* levels include antibiotics, cough medicines, laxatives, and steroids.
- ✚ Drugs that may cause *decreased* levels include diuretics (e.g., furosemide) and steroids.

### Procedure and patient care

- See inside front cover for Routine Urine Testing.
- If  $FE_{Na}$  is ordered, collect a venous blood sample in a gold-top tube for serum creatinine and sodium measurements.

### Abnormal findings

#### ▲ Increased levels

Dehydration  
Starvation  
Adrenocortical insufficiency  
Diuretic therapy  
Hypothyroidism  
SIADH secretion  
Diabetic ketoacidosis  
Toxemia of pregnancy

#### ▼ Decreased levels

Congestive heart failure  
Malabsorption  
Diarrhea  
Renal failure  
Cushing syndrome  
Aldosteronism  
Diaphoresis  
Pulmonary emphysema  
Inadequate sodium intake

notes

## spinal x-ray (Cervical, thoracic, lumbar, sacral, and coccygeal x-ray studies)

**Type of test** X-ray

### Normal findings

Normal spinal vertebrae

### Test explanation and related physiology

Spinal x-ray studies may be performed to evaluate any area of the spine. They usually include anteroposterior, lateral, and oblique views of these structures. These x-ray images are often done to assess back or neck pain, degenerative arthritic changes, traumatic fractures, tumor metastasis, spondylosis (stress fracture of the vertebrae), and spondylolisthesis (slipping of one vertebral disc on the other). Cervical spinal x-ray studies are performed in cases of multiple trauma to ensure that there is no fracture before the patient is moved or the neck is manipulated. However, CT scanning of the cervical vertebrae is increasingly becoming the standard of practice to ensure that there is no cervical fracture. Spinal x-rays are very helpful in evaluating children and adults for spinal alignment abnormalities (e.g., kyphosis, scoliosis). Magnetic resonance imaging is another very accurate method of evaluating the spine.

### Contraindications

- Patients who are pregnant unless the benefits outweigh the risks

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Instruct the patient to remove any metal objects covering the area to be visualized.
  - Immobilize the patient if a spinal fracture is suspected. Apply a neck brace if a cervical spine fracture is suspected.
- PT** Tell the patient that no fasting or sedation is required; however, if a fracture is suspected, the patient may be kept NPO.

#### During

- Note that the patient is placed on an x-ray table. Anterior, posterior, lateral, and oblique x-ray images are taken of the desired area on the spinal cord. These same views can also be obtained with the patient in the standing position.

## 840 spinal x-ray

- Note that a radiologic technologist takes spinal x-ray images in a few minutes.
- ▶ Tell the patient that no discomfort is associated with this study.

### After

- Note that positioning and patient activity depend on test results.

### Abnormal findings

Degenerative arthritis changes

Traumatic or pathologic fracture

Spondylosis

Spondylolisthesis

Metastatic tumor invasion

Scoliosis

Suspected spinal osteomyelitis

### notes

## sputum culture and sensitivity (C&S, Culture and Gram stain)

**Type of test** Sputum

### Normal findings

Normal upper respiratory tract

### Test explanation and related physiology

Sputum cultures are obtained to determine the presence of pathogenic bacteria in patients with respiratory infections (e.g., pneumonia). A *Gram stain* is the first step in the microbiologic analysis of sputum. Staining of sputum provides an opportunity to classify bacteria as gram positive or gram negative. This may be used to guide drug therapy until the culture and sensitivity (C&S) report is complete. The sputum sample is then applied to a series of bacterial culture plates. The bacteria that grow on those plates 1 to 3 days later are then identified. Determinations of bacterial sensitivity to various antibiotics are done to identify the most appropriate antimicrobial drug therapy. This is done by observing a ring of growth inhibition around an antibiotic plug in the culture medium.

Sputum for C&S should be collected before antimicrobial therapy is initiated unless the test is being performed to evaluate the effectiveness of medications already being given. Preliminary reports are usually available in 24 hours. Cultures require at least 48 hours for completion. Sputum cultures for fungus and *Mycobacterium tuberculosis* may take 6 to 8 weeks.

### Procedure and patient care

#### Before

- PT** Explain the procedure for sputum collection to the patient.
- PT** Remind the patient that sputum must be coughed up from the lungs and that saliva is not sputum.
  - Hold antibiotics until after the sputum has been collected.
  - If an elective specimen is to be obtained, give the patient a sterile sputum container on the night before the sputum is to be collected so that the morning specimen may be obtained on arising.
- PT** Instruct the patient to rinse out his or her mouth with water before the sputum collection to decrease contamination of the sputum by particles in the oropharynx.

## 842 sputum culture and sensitivity

### During

- Note that sputum specimens are best when the patient first awakens in the morning and before eating or drinking.
- Collect at least 1 teaspoon of sputum in a sterile sputum container.
- Usually obtain sputum by having the patient cough after taking several deep breaths.
- If the patient is unable to produce a sputum specimen, stimulate coughing by lowering the head of the patient's bed or giving the patient an aerosol administration of a warm hypertonic solution.
- Note that other methods to collect sputum include endotracheal aspiration, fiberoptic bronchoscopy, and transtracheal aspiration.

### After

- PT** Inform the patient to notify the nurse as soon as the sputum is collected.
- Label the sputum and send it to the laboratory as soon as possible.
- Note any current antibiotic therapy on the laboratory slip.

### Abnormal findings

Bacterial infection (e.g., pneumonia)

Viral infection

Atypical bacterial infection (e.g., tuberculosis)

### notes

## sputum cytology

**Type of test** Sputum

### Normal findings

Normal epithelial cells

### Test explanation and related physiology

Tumors in the pulmonary system frequently slough cells into the sputum. When the sputum is gathered, the cells are examined. If the cytologic test result is positive, malignant cells are seen, indicating a lung tumor. If only normal epithelial cells are seen, either no malignancy exists or any existing tumor is not shedding cells. Therefore a positive test result indicates malignancy; a negative test result means nothing.

Bronchoscopy and percutaneous lung biopsy have supplanted the need for sputum cytology to a large degree. Now its greatest use is in patients who have an abnormal chest x-ray result, productive cough, and nothing visible on bronchoscopy.

### Procedure and patient care

#### Before

- PT** Explain the procedure for sputum collection to the patient.
- PT** Remind the patient that sputum must be coughed up from the lungs and that saliva is not sputum.
- Give the patient a sterile sputum container on the night before the sputum is to be collected so that the morning specimen may be obtained on arising.

#### During

- Sputum specimens are collected as described on p. 841.
- Usually collect sputum for cytology on three separate occasions.

#### After

- PT** Instruct the patient to notify the nurse as soon as the sputum is collected.
- Label the specimen, and send it to the laboratory as soon as possible.

### Abnormal findings

Malignancies

notes

**stool cancer screening** (stool for occult blood, Stool for OB, Fecal occult blood test [FOBT], Fecal immunotest [FIT], DNA stool sample)

**Type of test** Stool

### Normal findings

No occult blood within stool

### Test explanation and related physiology

This test is used for colorectal cancer screening of asymptomatic individuals. Normally only minimal quantities (2-2.5 mL) of blood are passed into the GI tract. Usually this bleeding is not significant enough to cause a positive result in stool for occult blood (OB) testing. This test can detect OB when as little as 5 mL of blood is lost per day.

Tumors of the intestine grow into the lumen and are subjected to repeated trauma by the fecal stream. Eventually the friable neovascular tumor ulcerates and bleeding occurs. Most often, bleeding is so slight that gross blood is not seen in the stool. The blood can be detected by chemical assay or by immunohistochemistry. Guaiac is the most commonly performed chemical assay.

OB can also be detected by immunochemical methods that detect the human globin portion of hemoglobin using monoclonal antibodies. These tests are called *fecal immunochemical test (FIT)* or *immunochemical fecal occult blood test (iFOBT)*. These methods are as sensitive as guaiac testing but are not affected by red meats or plant oxidizers as described in the Interfering factors section. Immunochemical methods may fail to recognize OB from the upper GI tract because the globin is digested by the time it gets in the stool.

The *DNA stool sample test* is more sensitive than guaiac testing in the detection of significant colorectal precancerous, benign, and malignant tumors. Because most precancerous polyps do not bleed, they can be missed by FOBT. In contrast, all precancerous polyps shed cells that contain abnormal DNA. So, a stool-based DNA test designed to detect this DNA promises to be more accurate in the detection of precancerous polyps, which, when detected, can be removed before they turn into cancer. The test is an easy-to-use home kit to collect a stool sample and mail it to a laboratory for analysis. The test checks for DNA changes that could indicate cancer or precancerous polyps and also checks for the presence of blood in the stool that can indicate cancer.

Benign and malignant GI tumors, ulcers, inflammatory bowel disease, arteriovenous malformations, diverticulosis, and



hematobilia (hemobilia) can all cause OB in the stool. Other more common abnormalities (e.g., hemorrhoids, swallowed blood from oral or nasopharyngeal bleeding) may also cause OB in the stool.

When OB testing is properly performed, a positive result obtained on multiple specimens collected on successive days warrants a thorough GI evaluation, usually with esophagogastroduodenoscopy (EGD) (see p. 382) and colonoscopy (see p. 258). Regular screening, beginning at age 50 years, can reduce the number of people who die of colorectal cancer by as much as 60%.

Reducing or oxidizing agents (e.g., iron, radish, cantaloupe, cauliflower, vitamin C) can affect the results of guaiac or FIT. Furthermore, neither FIT nor guaiac testing detects slow upper GI bleeding because globin and heme are degraded during intestinal transit.

### Interfering factors

- Vigorous exercise
- Bleeding gums after a dental procedure
- Ingestion of red meat within 3 days before testing
- Ingestion of peroxidase-rich fruits and vegetables (turnips, artichokes, mushrooms, radishes, horseradishes, broccoli, bean sprouts, cauliflower, oranges, bananas, cantaloupes, and grapes) may affect results.
- Diarrhea or blood in the urine or stool (e.g., bleeding hemorrhoids, bleeding cuts or wounds on their hands, rectal bleeding, or menstruation)
- ✎ Drugs that may cause GI bleeding include anticoagulants, aspirin, colchicine, iron preparations (large doses), nonsteroidal antiarthritics, and steroids.
- ✎ Drugs that may cause false-positive results include colchicine, iron, oxidizing drugs (e.g., iodine, bromides, boric acid), and rauwolfia derivatives.
- ✎ Drugs that may cause false-negative results include vitamin C.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to refrain from eating any red meat for at least 3 days before the test.
- PT** Instruct the patient to refrain from drugs known to interfere with OB testing.
- PT** Instruct the patient as to the method of obtaining appropriate stool specimens. Tests may be done at home with specimen cards (Hemoccult) and mailed.

## 846 stool cancer screening

- PT** Instruct the patient not to mix urine with the stool specimen.
- PT** Inform the patient as to the need for multiple specimens obtained on separate days to increase the test's accuracy.
  - Note that in some lefts a high-residue diet is recommended to increase the abrasive effect of the stool.
  - Be gentle in obtaining stool by digital rectal examination. A traumatic digital examination can cause a false-positive stool, especially in patients with prior anorectal disease, such as hemorrhoids.

### During

#### *Hemoccult slide test*

- Place a stool sample on one side of guaiac paper.
- Place two drops of developer on the other side.
- Note that bluish discoloration indicates OB in the stool.

#### *Tablet test*

- Place a stool sample on the developer paper.
- Place a tablet on top of the stool specimen.
- Put two or three drops of tap water on the tablet and allow it to flow onto the paper.
- Note that bluish discoloration indicates OB in the stool.

#### *DNA Home Test*

- Place the bracket on the toilet.
- Add container to the bracket.
- Have bowel movement.
- Place a small stool sample in the smaller tube.
- Place preservative on the stool in the larger container of stool.
- Replace the top on the container and mail both the container and the smaller tube to the address on the enclosed label.

### After

- PT** Inform the patient of the results.
  - If the tests are positive, inquire whether the patient violated any of the preparation recommendations.

### Abnormal findings

GI tumor	Ischemic bowel disease
Polyps	GI trauma
Ulcer	Recent GI surgery
Varices	Hemorrhoids
Inflammatory bowel disease	Esophagitis
Diverticulosis	Gastritis

### notes

## stool culture (Stool for culture and sensitivity [C&S], Stool for ova and parasites [O&P])

**Type of test** Stool

### Normal findings

Normal intestinal flora

### Test explanation and related physiology

Normally, stool contains many bacteria and fungi. The more common organisms include *Enterococcus* spp., *Escherichia coli*, *Proteus* spp., *Pseudomonas* spp., *Staphylococcus aureus*, *Candida albicans*, *Bacteroides* spp., and *Clostridium* spp. Bacteria are indigenous to the bowel. Sometimes normal stool flora can become pathogenic if overgrowth of the bacteria occurs as a result of antibiotics (e.g., *Clostridium difficile*), immunosuppression, or overaggressive catharsis. *Salmonella* spp., *Shigella* spp., *Campylobacter* spp., *Yersinia* spp., pathogenic *E. coli*, *Clostridium*, and *Staphylococcus* spp. are acquired bacteria that can infect the bowel. Parasites also may affect the stool. Common parasites are *Ascaris* spp. (hookworm), *Strongyloides* spp. (tapeworm), and *Giardia* spp. (protozoans). Identification of any of these pathogens in the stool incriminates that parasite as the etiology of the infectious enteritis.

Infections of the bowel from bacteria, virus, or parasites usually present as diarrhea, excessive flatus, and abdominal discomfort. Patients who have been drinking well water, have been on prolonged antibiotics, or have traveled outside of the United States are especially susceptible.

### Interfering factors

- Urine may inhibit the growth of bacteria. Therefore urine should not be mixed with the feces during collection of a stool sample.
- Recent barium studies may obscure the detection of parasites.
- ✎ Drugs that may affect test results include antibiotics, bismuth, and mineral oil.

### Procedure and patient care

#### Before

- PT** Explain the method of stool collection to the patient. Be matter of fact to avoid any embarrassment to the patient.

## 848 stool culture

- PT** Instruct the patient not to mix urine or toilet paper with the stool specimen.
- PT** Instruct the patient to use an appropriate collection container.

### During

- PT** Instruct the patient to defecate into a clean bedpan.
  - Place a small amount of stool in a sterile collection container.
  - Send mucus and blood streaks with the specimen.
  - If a rectal swab is to be used, wear gloves and insert the cotton-tipped swab at least 1 inch into the anal canal. Then rotate the swab for 30 seconds and place it into the clean container.

### Tape test

- Use this test when pinworms (*Enterobius*) are suspected.
- Place clear tape in the patient's perianal region. (This is especially helpful in children.)
- Because the female worm lays her eggs at night around the perianal area, apply the tape before bedtime and remove it in the morning before the patient gets out of bed.
- Press the sticky surface of the tape directly to a glass slide and examine microscopically for pinworm ova.

### After

- Handle the stool specimen carefully as though it were capable of causing infection.
- Indicate on the laboratory slip any antibiotics that the patient may be taking.
- Promptly send the stool specimen to the laboratory. Delays in transfer of the specimen may affect viability of the organism.
- Note that some enteric pathogens occasionally take as long as 6 weeks to isolate.
- When pathogens are detected, maintain isolation of the patient's stool until therapy is completed.

### Abnormal findings

Bacterial enterocolitis  
Protozoan enterocolitis  
Parasitic enterocolitis

### notes

## stool for leukocytes (White blood cell stool test, Fecal leukocyte stain, Stool for white cells)

**Type of test** Stool

### Normal findings

≤ 2/hpf

### Test explanation and related physiology

Leukocytes are not normally seen in stools in the absence of infection or inflammatory bowel diseases. Fecal leukocytosis is a response to infection with microorganisms that invade tissue or produce toxins. This leads to tissue damage of the bowel wall, generating a vigorous leukocyte infiltration. This test is used to identify intestinal infections, diarrheal diseases, or inflammatory bowel diseases.

Fecal leukocytes are commonly found in patients with bacterial infections such as *Shigella*, *Campylobacter*, *Salmonella*, *Yersinia*, or *Clostridium* spp. Amebiasis is also associated with fecal leukocytes. The greater the number of leukocytes, the greater the likelihood of infection. Diarrhea caused by most parasites (e.g., *Giardia* spp.) or by viral infections (e.g., Norwalk, rotaviruses, or adenoviruses) do not cause leukocytes in the stool. Therefore a negative fecal leukocyte test result does not rule out other potential problems.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting is needed.

#### During

- Collect a random stool specimen at least the size of a walnut. Liquid stool may be used.
- Carefully follow instructions on container. Fresh, ECOFIX-preserved, or polyvinyl alcohol-preserved stool must be sent to the laboratory.

#### After

- PT** Tell the patient that results will be available within 1 to 2 days.

### Abnormal findings

Infectious colitis such as cholera, *C. difficile*, *Salmonella* spp., and so on.

*Inflammatory Bowel diseases such as* ulcerative colitis, or Crohn disease

***Streptococcus* serologic testing** (Antistreptolysin O titer [ASO], Antideoxyribonuclease-B titer, [Anti-DNase-B, ADNase-B, ADB], *Streptococcus* group B antigen detection, Streptozyme)

**Type of test** Blood, cerebrospinal fluid (CSF)

**Normal findings**

Antistreptolysin O titer

Adult/elderly:  $\leq 160$  Todd units/mL

Child:

Newborn: similar to mother's value

6 months-2 years:  $\leq 50$  Todd units/mL

2-4 years:  $\leq 160$  Todd units/mL

5-12 years: 170-330 Todd units/mL

Antideoxyribonuclease-B titer

Adult:  $\leq 85$  Todd units/mL or titer  $\leq 1:85$

Child:

Preschool age:  $\leq 60$  Todd units/mL or titer  $\leq 1:60$

School age:  $\leq 170$  Todd units/mL or titer  $\leq 1:170$

Streptozyme

Titer  $< 1:100$

*Streptococcus* group B antigen

None detected

**Test explanation and related physiology**

Infection by group A *Streptococcus* is unique because it can be followed by a serious nonpurulent complication (e.g., rheumatic fever, scarlet fever, glomerulonephritis). Serologic tests are used primarily to determine whether a previous group A *Streptococcus* infection (pharyngitis, pyoderma, pneumonia) has caused a poststreptococcal disease. These poststreptococcal diseases occur after the infection and after a period of latency during which the patient is asymptomatic.

These antibodies are directed against streptococcal extracellular products that are primarily enzymatic proteins. Serial rising titers of these antibodies over several weeks followed by a slow fall in titers are more supportive of the diagnosis of a previous streptococcal infection than is a single titer.

One such extracellular enzyme produced by streptococcus is called *streptolysin O*, which has the ability to destroy (lyse) RBC corpuscles. The streptolysin O is antigenic, stimulating the immunologic production of a neutralizing ASO antibody. ASO appears in the serum 1 week to 1 month after the onset of a streptococcal infection. A high ASO titer is not specific for

a certain type of poststreptococcal disease (i.e., rheumatic fever versus glomerulonephritis) but merely indicates that a streptococcal infection is or has been present.

Like the ASO titer, *ADB* is used to detect previous streptococcal infections. Although this test may be more sensitive than the ASO titer, it is not used alone in the evaluation of streptococcal infections because its results are too variable.

The *Streptozyme* assay detects antibodies to multiple extracellular antigens of group A *Streptococcus*, including antistreptolysin O, antistreptokinase, and antihyaluronidase. Approximately 80% of specimens positive by streptozyme have antistreptolysin O, and 10% have antistreptokinase and/or antihyaluronidase. The remaining 10% of positive samples are apparently caused by *ADB* antibodies or other streptococcal extracellular antigens.

*Streptococcus group B antigens* accumulate in CSF, serum, or urine and provide a direct qualitative detection of bacterial antigens. These antigens indicate acute infection and are not related to poststreptococcal sequelae as described previously. Confirmatory diagnosis of streptococcal infection is done by cultures (see p. 872).

### Interfering factors

- Increased beta-lipoprotein levels inhibit streptolysin O and give a falsely high ASO titer.
- ✚ Drugs that may cause *decreased* ASO levels include adrenocorticosteroids and antibiotics.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

#### ▲ Increased levels

Streptococcal infection  
 Acute rheumatic fever  
 Acute glomerulonephritis  
 Bacterial endocarditis  
 Scarlet fever  
 Streptococcal pyoderma

### notes

**substance abuse testing** (Urine drug testing, Drug screening, Toxicology screening)

**Type of test** Urine; blood; various

**Normal findings**

Negative

**Test explanation and related physiology**

Substance abuse testing is used mostly by employers and law enforcement agencies. Industrial testing is used at the time of preemployment, prepromotion, annual physical examination, after an accident when there is reasonable suspicion, or for random testing or follow-up treatment surveillance.

Most commonly, a drug screen is performed to detect small amounts of any number of metabolites of commonly used drugs. If the screen result is positive, a more accurate and quantitative test is performed on the same specimen. Drug screens are available for a variety of drug categories. The most common are amphetamines, barbiturates, benzodiazepines, cannabinoids (marijuana [THC]), cocaine, methamphetamine, opiates (morphine and heroin), phencyclidine (PCP), carisoprodol, meprobamate, and propoxyphene (Table 36). Alcohol testing is most commonly used by law enforcement (see ethanol, p. 391). Not only is drug testing helpful in identifying users, it also acts as a deterrent. Athletes can be tested for anabolic hormones, stimulants, diuretics, beta-blockers, street drugs, antiestrogens, erythropoietin, and beta-2 agonists that may unfairly improve their performance. Health and life insurance companies routinely test for illicit drugs.

Until recently, substance abuse testing has used urine exclusively as the sample of choice. Urine drug testing is generally inexpensive. Urine is easily obtained, and it contains a large amount of drug metabolites. More important, whereas urine can identify drug usage for several days after the last usage, blood testing reflects drug usage only during the past few hours.

Saliva, breath, hair, and sweat are becoming increasingly important and accurate specimens for specific drug testing. These testing methods are very expensive, however. Hair samples detect the presence of drugs used during the past 3 months. In addition, hair and nail samples may be used to detect or document exposure to arsenic and mercury.



**TABLE 36** Typical multipanel drug screen

Drugs or drug classes	Screen	Confirmation <sup>a</sup>
Marijuana	20 ng/mL	5 ng/mL
Cocaine	150 ng/mL	50 ng/mL
Opiates	300 ng/mL	5 ng/mL
Oxycodone	100 ng/mL	5 ng/mL
Phencyclidine	25 ng/mL	10 ng/mL
Amphetamines	300 ng/mL	200 ng/mL
MDMA (Ecstasy)	500 ng/mL	200 ng/mL
Barbiturates	200 ng/mL	50 ng/mL
Benzodiazepines	200 ng/mL	20 ng/mL
Methadone	150 ng/mL	10 ng/mL
Propoxyphene	300 ng/mL	10 ng/mL

<sup>a</sup> Confirmatory tests are more sensitive and can detect metabolites at lower levels.

Because a positive result can have a profound effect on a person’s life, job, and accountability, it is common for a drug abuser to attempt to alter the urine specimen. Therefore the urine sample is tested for odor, color, temperature, creatinine, pH, and specific gravity to ensure that it is a proper specimen.

*Toxicology screening* tests for drug overdose and poisoning (e.g., lead and carbon monoxide) are best performed on blood. Results indicate current drug levels, which are used to determine or alter therapy. Toxicology studies are used to incriminate drugs as a cause or factor in the death of a person.

### Interfering factors

- Poppy seeds can cause positive opiate results.
- Secondhand marijuana smoke can cause positive THC results.
- Detergents, bicarbonates, salt tablets, and blood can all foil accurate drug testing in a urine specimen.
- ✚ Ibuprofen can cause a false-positive THC result in some assays.
- ✚ Cold remedies can cause false-positive amphetamine results in some assay systems.
- ✚ Antibiotics (e.g., amoxicillin) can cause false-positive results for heroin and/or cocaine.
- ✚ Aggressive use of diuretics can *decrease* urine drug levels.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient or significant others.
- If the specimen is obtained for medicolegal testing, ensure that the patient or family member has signed a consent form.
  - Obtain a list of prescription medicines that the patient is taking that may alter or confuse screening results. Obtain as much information as possible about the drug type, amount, and ingestion time.
  - Carefully assess the patient for respiratory distress.

### During

- Collect blood and urine samples as designated by the laboratory.
- Ensure that patients provide their own urine. Urine specimens for substance abuse testing are usually collected in the presence of a healthcare provider.
- Be sure that the patient does not alter the urine specimen.
- For hair testing, cut 50 strands of hair from the scalp.
- A second confirmatory specimen may be obtained (and is used if results are positive).
- Collect gastric contents as indicated by the specific institution. A nasogastric tube is required.

### After

- Apply pressure to the venipuncture site.
- Refer the patient for appropriate drug and psychiatric counseling.
- Follow the chain of custody for the specimen as provided by standard guidelines of the institution.
- Place the specimen in the required container for delivery.
- Check the temperature of urine specimens within 3 minutes after voiding. The temperature should be between 97° and 99° F.
- The specimen may be sent to a nationally certified laboratory for federal workers or workplace testing. Local hospital laboratories are often able to test for many drugs.

## Abnormal findings

Positive drug level

notes

## swallowing examination (Videofluoroscopy swallowing examination)

**Type of test** X-ray with contrast dye

### Normal findings

Normal swallowing function and complete clearing of radiographic material through the upper digestive tract

### Test explanation and related physiology

This test is performed to identify problems that exist in a patient who is unable to swallow. Problems in swallowing may result from local structural diseases such as tumors, upper esophageal diverticula, inflammation, extrinsic compression of the upper GI tract, or surgery on the oropharyngeal tract. Motility disorders of the upper GI tract (e.g., Zenker diverticulum) and neurologic disorders (e.g., stroke syndrome), Parkinson disease, and neuropathies also may cause difficulty in swallowing. Videofluoroscopy of the swallowing function allows a speech pathologist to delineate more clearly the exact pathology in the swallowing mechanism. This procedure then can be used to determine the most appropriate treatment and teach the patient the proper swallowing technique.

This test is performed by asking the patient to swallow barium or a barium-containing meal. With the use of videofluoroscopy, the swallowing function is visualized and documented. Morphologic abnormalities and functional impairment can be identified easily using the slow-frame progression and reversal that is available with videofluoroscopy. Although this test is similar to the barium swallow (see p. 131), finer details of swallowing can be evaluated with the use of videofluoroscopy.

### Contraindications

- Patients who aspirate their saliva are not candidates for this swallowing examination because they will require nonswallowing methods of alimentation.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Explain to the patient that no preparation is required.

## 856 swallowing examination

### During

- In the radiology department, the patient is asked to swallow a barium-containing meal. The consistency of the meal will be determined by the speech therapist and radiologist. The food may be liquid, semisoft (e.g., applesauce), or solid (e.g., a tea biscuit). While the patient is swallowing, videofluoroscopy is recorded in both the lateral and the anterior positions.
- The video is then repeatedly examined and reexamined by the radiologist and speech pathologist.

### After

- No catharsis is required.

### Abnormal findings

Oral pharyngeal inflammation

Cancer

Extrinsic compression

Neuromuscular disorder

Achalasia

Upper GI motility disorder (e.g., stroke syndrome, Parkinson disease, peripheral neuropathy)

Diffuse esophageal spasms

Zenker diverticulum

notes

## sweat electrolytes test (Iontophoretic sweat test)

**Type of test** Fluid analysis

### Normal findings

Sodium values in children

Normal: < 70 mEq/L

Abnormal: > 90 mEq/L

Equivocal: 70-90 mEq/L

Chloride values in children

Normal: < 50 mEq/L

Abnormal: > 60 mEq/L

Equivocal: 50-60 mEq/L

### Test explanation and related physiology

Patients with cystic fibrosis (CF) have increased sodium and chloride contents in their sweat. This forms the basis of this test, which is both sensitive and specific for CF. CF is an inherited disease characterized by abnormal secretion by exocrine glands in the bronchi, small intestines, pancreatic ducts, bile ducts, and skin (sweat glands). Sweat induced by electrical current (*pilocarpine iontophoresis*) is collected, and its sodium and chloride contents are measured. The degree of abnormality is no indication of the severity of CF; it merely indicates that the patient has the disease.

This test is also used to screen children or siblings of CF patients for the disease. Almost all patients with CF have sweat sodium and chloride contents two to five times greater than normal values. In patients with suspicious clinical manifestations, these levels are diagnostic of CF.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient and parents.

**PT** Tell the patient and parents that no fasting is required.

#### During

- Note the following procedural steps:
  1. For iontophoresis, a low-level electrical current is applied to the test area (the thigh in infants and the forearm in older children).
  2. The positive electrode is covered by gauze and saturated with pilocarpine hydrochloride, a stimulating drug that induces sweating.

## 858 sweat electrolytes test

3. The negative electrode is covered by gauze saturated with a bicarbonate solution.
  4. The electrical current is allowed to flow for 5 to 12 minutes.
  5. The electrodes are removed, and the arm is washed with distilled water.
  6. Paper disks are placed over the test site with the use of clean, dry forceps.
  7. These disks are covered with paraffin to obtain an airtight seal, preventing evaporation of sweat.
  8. After 1 hour the paraffin is removed. The paper disks are transferred immediately by forceps to a weighing jar and sent for sodium and chloride analysis.
  9. A *screening test* may be done to detect sweat chloride levels. For screening, a test paper containing silver nitrate is pressed against the child's hand for several seconds. The test is positive when the excess chloride combines with the silver nitrate to form white-silver chloride on the paper (i.e., the child with CF will leave a "heavy" handprint on the paper).
  10. A positive screening test is usually validated by iontophoresis.
- Note that an experienced technologist performs the sweat test in approximately 90 minutes in the laboratory or at the patient's bedside.
- PT** Inform the patient that the electrical current is small and no discomfort or pain is generally associated with this test.

### After

- PT** Initiate extensive education and counseling for the patient and parents if the results indicate CF.

### Abnormal findings

CF

notes

**syphilis detection test** (Serologic test for syphilis [STS], Venereal Disease Research Laboratory [VDRL], Rapid plasma reagin [RPR], Fluorescent treponemal antibody test [FTA])

**Type of test** Blood

### Normal findings

Negative or nonreactive

### Test explanation and related physiology

Serologic tests are used to diagnose and to document successful therapy for syphilis (see p. 813). Syphilis is caused by the spirochete *Treponema pallidum* that cannot be isolated in culture. Two groups of antibodies form the basis for these tests. The first and older of these tests detects a nontreponemal antibody called *reagin*, which reacts to phospholipids similar to lipids in the membrane of *T. pallidum*. The nontreponemal antibody tests are relatively nonspecific and lack sensitivity. These antibodies would be detected by the *Wassermann test*, *VDRL test*, or *RPR test*. These test results become positive after 2 weeks from the patient's inoculation with *T. pallidum* and return to normal after adequate treatment is administered. The test result is positive in nearly all primary and secondary stages of syphilis and in two-thirds of patients with tertiary syphilis. Screening for syphilis is usually done during the first prenatal checkup for pregnant women using the VDRL or RPR. VDRL is the only test that can be used on cerebral spinal fluid (CSF) when evaluating neurosyphilis. They are also used to document the success of treatment.

If these nontreponemal serologic test results are positive, the diagnosis must be confirmed by the second type of syphilis test, called *Treponema test*, such as the *FTA absorption test (FTA-ABS)* or the *microhemagglutination assay (MHA-TP)*. These tests for a more specific antibody are more accurate than the VDRL and RPR tests. The FTA-ABS and MHA-TP are technically simple to perform, but they are labor intensive and require subjective interpretation by testing personnel. In contrast, the syphilis IgG *enzyme immunoassay (EIA)* is a treponemal test for the detection of IgG class antibodies.

During early primary syphilis, the first antibodies to appear are IgM, with IgG antibodies reaching significant titers later in the primary phase. As the disease progresses into the secondary phase, IgG *T. pallidum* antibodies reach peak titers. *T. pallidum* IgG antibodies persist indefinitely, regardless of the course of the disease. If syphilis IgG and/or IgM is positive, results can

## 860 syphilis detection test

be confirmed with FTA or MHA testing. The IgG- and IgM-specific antibodies assist in determining the etiology of neonatal syphilis. IgM does not pass through the placenta and, if positive, indicates active neonatal infection.

### Interfering factors

- Excessive hemolysis and gross lipemia may affect test results.
- Excess chyle in the blood may interfere with the test results.
- Many conditions cause false-positive results when VDRL and RPR tests are used. Some of these conditions include *Mycoplasma pneumoniae*, malaria, acute bacterial and viral infections, autoimmune diseases, and pregnancy.
- Recent ingestion of alcohol may alter the test results.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
  - Fasting: verify with laboratory
  - Blood tube commonly used: red
  - Check with the laboratory regarding fasting requirements. Some prefer collecting the specimen before meals. Some laboratories request that the patient refrain from alcohol for 24 hours before the blood test.
- PT** the test result is positive, instruct the patient to inform recent sexual contacts so that they can be evaluated.
- PT** If the test result is positive, be sure the patient receives the appropriate antibiotic therapy.

### Abnormal findings

Syphilis

notes



**testosterone** (Dihydrotestosterone [DHT])**Type of test** Blood**Normal findings****Free testosterone, pg/mL****Male**

Tanner stage I:  $\leq 3.7$   
 Tanner stage II: 0.3-21  
 Tanner stage III: 1.0-98  
 Tanner stage IV: 35.0-169  
 Tanner stage V: 41.0-239

**Female**

Postmenopausal: 0.6-3.8  
 Tanner stage I:  $< 2.2$   
 Tanner stage II: 0.4-4.5  
 Tanner stage III: 1.3-7.5  
 Tanner stage IV: 1.1-15.5  
 Tanner stage V: 0.8-9.2

**% Free testosterone**

Adult male: 1.6%-2.9%  
 Adult female: 0.1%-0.3%

**Total testosterone, ng/dL****Age****Male****Female**

7 months-9 years (Tanner stage I)	$< 30$	$< 30$
10-13 years (Tanner stage II)	$< 300$	$< 40$
14-15 years (Tanner stage III)	170-540	$< 60$
16-19 years (Tanner stage IV, V)	250-910	$< 70$
20 years and above	280-1080	$< 70$

**Dihydrotestosterone**

Adult male: 240-650 pg/mL  
 Adult female:  $\leq 300$  pg/mL

**Test explanation and related physiology**

Testosterone levels are used to evaluate ambiguous sex characteristics, precocious puberty, virilizing syndromes in the female, and infertility in the male. This test can also be used as a tumor marker for rare tumors of the ovary and testicle.

Androgens include dehydroepiandrosterone (DHEA), androstenedione, and testosterone. DHEA is produced in the adrenal glands during cortisol and aldosterone formation, and it is also produced *de novo* by the testes or the ovaries. DHEA is the precursor of androstenedione, which is the precursor of testosterone (and estrogen).

Testosterone levels vary by sex and stage of maturity (indicated by Tanner stage). In men, most of the testosterone is made by the Leydig cells in the testicle; this accounts for 95% of the circulating testosterone in men. In women, about half of the testosterone is made by the conversion of DHEA to testosterone in the peripheral fat tissue. Another 30% is made by the same conversion of DHEA in the adrenal gland, and 20% is made directly by the ovaries.

Physiologically, testosterone stimulates spermatogenesis and influences the development of male secondary sex characteristics. Overproduction of this hormone in a young male may cause precocious puberty. This can be caused by testicular, adrenal, or pituitary tumors. Overproduction of this hormone in females causes masculinization, which is manifested as amenorrhea and excessive growth of body hair (hirsutism). Ovarian and adrenal tumors/hyperplasia and medications (e.g., danazol) are all potential causes of masculinization in the female. Reduced levels of testosterone in the male suggest hypogonadism or Klinefelter syndrome.

Dihydrotestosterone (DHT) is the principal androgen made in body tissues, particularly the prostate. Levels of DHT remain normal with aging, despite a decrease in the plasma testosterone, and are not elevated in benign prostatic hyperplasia. Measurement of this hormone is useful in monitoring patients receiving 5 alpha-reductase inhibitor therapy, such as finasteride or chemotherapy, which may affect prostate function. It is also useful in evaluating patients with possible 5 alpha-reductase deficiency.

There are several *testosterone stimulation tests* that can be performed to more accurately evaluate hypogonadism. Human chorionic gonadotropin, clomiphene, and GnRH can be used to stimulate testosterone secretion.

*17-Ketosteroids (17-KS)* are metabolites of the testosterone and nontestosterone androgenic sex hormones that are excreted in the urine.

### Interfering factors

- ✦ Drugs that may cause *increased* testosterone levels include anticonvulsants, barbiturates, estrogens, and oral contraceptives.
- ✦ Drugs that may cause *decreased* testosterone levels include alcohol, androgens, dexamethasone, diethylstilbestrol, digoxin, ketoconazole, phenothiazine, spironolactone, and steroids.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Because testosterone levels are the highest in the early morning hours, blood should be drawn in the morning.

## Abnormal findings

### ▲ Increased levels (male)

Idiopathic sexual precocity  
 Pinealoma  
 Encephalitis  
 Congenital adrenal hyperplasia  
 Adrenocortical tumor  
 Testicular or extragonadal tumor  
 Hyperthyroidism  
 Testosterone resistance syndromes

### ▼ Decreased levels (male)

Klinefelter syndrome  
 Cryptorchidism  
 Primary and secondary hypogonadism  
 Trisomy 21 (Down syndrome)  
 Orchidectomy  
 Hepatic cirrhosis

### ▲ Increased levels (female)

Ovarian tumor  
 Adrenal tumor  
 Congenital adrenocortical hyperplasia  
 Trophoblastic tumor  
 Polycystic ovaries  
 Idiopathic hirsutism

notes

## thoracentesis and pleural fluid analysis (Pleural tap)

**Type of test** Fluid analysis

### Normal findings

Gross appearance: Clear, serous, light yellow, 50 mL

Red blood cells (RBCs): None

White blood cells (WBCs): < 300/mL

Protein: < 4.1 g/dL

Glucose: 70-100 mg/dL

Amylase: 138-404 units/L

Alkaline phosphatase

Adult male: 90-240 units/L

Female: < 45 years: 76-196 units/L

Female: > 45 years: 87-250 units/L

Lactic dehydrogenase (LDH): Similar to serum LDH

Cytology: No malignant cells

Bacteria: None

Fungi: None

Carcinoembryonic antigen (CEA): < 5 ng/mL

### Test explanation and related physiology


Thoracentesis is an invasive procedure that entails insertion of a needle into the pleural space for removal of fluid ([Figure 40](#)). Pleural fluid is removed for diagnostic and therapeutic purposes. *Therapeutically*, it is done to relieve pain, dyspnea, and other symptoms of pleural pressure. Removal of this fluid also permits better radiographic visualization of the lung. *Diagnostically*, thoracentesis is performed to obtain and analyze fluid to determine the etiology of the pleural effusion.

Pleural fluid is classified according to transudate or exudate. This is an important differentiation and is very helpful in determining the etiology of the effusion. *Transudates* are most frequently caused by congestive heart failure, cirrhosis, nephrotic syndrome, and hypoproteinemia. *Exudates* are most often found in inflammatory, infectious, or neoplastic conditions. However, collagen vascular disease, pulmonary infarction, trauma, and drug hypersensitivity also may cause an exudative effusion.

Pleural fluid is usually evaluated for the following features.

#### Gross appearance

The color, optical density, and viscosity are noted as the pleural fluid appears in the aspirating syringe. Empyema is characterized by the presence of a foul odor and thick, puslike fluid. An

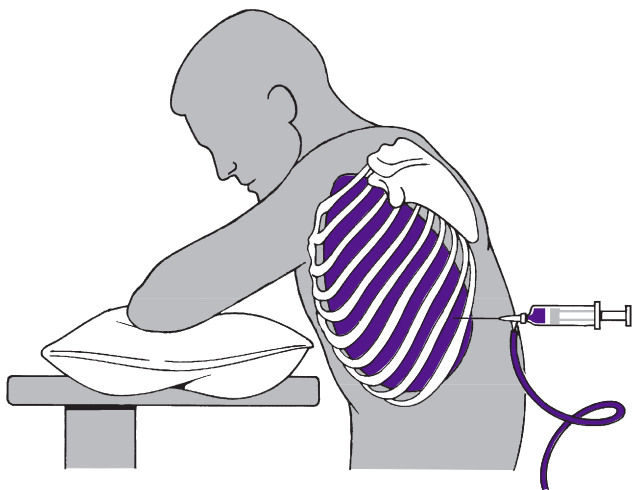

**FIG. 40** Thoracentesis. A needle is placed through the chest wall and into the fluid contained in the pleural cavity. A special one-way valve system is placed between the needle and the syringe to allow aspiration of fluid when the plunger of the syringe is pulled back and diversion of the fluid to a container when the plunger is pushed in.

opalescent, pearly fluid is characteristic of chylothorax (chyle in the pleural cavity).

### Cell counts

The WBC and differential counts are determined. A WBC count exceeding 1000/mL is suggestive of an exudate. The predominance of polymorphonuclear leukocytes usually is an indication of an acute inflammatory condition (e.g., pneumonia, pulmonary infarction, early tuberculosis [TB] effusion). When more than 50% of the WBCs are small lymphocytes, the effusion is usually caused by TB, or tumor. Normally, no RBCs should be present. The presence of RBCs may indicate neoplasms, TB, or intrathoracic bleeding.

### Protein content

Whereas total protein levels greater than 3 g/dL are characteristic of exudates, transudates usually have a protein content of less than 3 g/dL. The *albumin gradient* between serum and pleural fluid can differentiate better between the transudate and exudate

## 866 thoracentesis and pleural fluid analysis

natures of pleural fluid than can the total protein content. This gradient is obtained by subtracting the pleural albumin value from the serum albumin value. Values of 1.1 g/dL or more suggest a transudate. Values less than 1.1 g/dL suggest an exudate but will not differentiate the potential cause of the exudate (e.g., malignancy from infection or inflammation). The total protein ratio (fluid/serum) has been considered to be another accurate criterion differentiating transudate from exudate. A total protein ratio of fluid to serum of greater than 0.5 is considered to be an exudate.

### Lactic dehydrogenase

A pleural fluid/serum LDH ratio greater than 0.6 is typical of an exudate. An exudate is identified with a high degree of accuracy if the pleural fluid/serum protein ratio is greater than 0.5 and the pleural fluid/serum LDH ratio is greater than 0.6.

### Glucose

Usually pleural glucose levels approximate serum levels. Low values appear to be a combination of glycolysis by the extra cells and impairment of glucose diffusion because of damage to the pleural membrane. Values less than 60 mg/dL are occasionally seen in TB, or malignancy and typically occur in rheumatoid arthritis and empyema.

### Amylase

In a malignant effusion, the amylase concentration is slightly elevated. Amylase levels above the normal range for serum or two times the serum level are seen when the effusion is caused by pancreatitis or rupture of the esophagus associated with leakage of salivary amylase.

### Triglyceride

Measurement of triglyceride levels is an important part of identifying chylous effusions. These effusions are usually produced by obstruction or transection of the lymphatic system caused by lymphoma, neoplasm, trauma, or recent surgery. The triglyceride value in a chylous effusion exceeds 110 mg/dL.

### Gram stain and bacteriologic culture

These tests are routinely performed when bacterial pneumonia or empyema is a possible cause of the effusion. If possible, these should be done before initiation of antibiotic therapy.

### Cultures for *Mycobacterium tuberculosis* and fungus

Tuberculosis is less often a cause for pleural effusion in the United States today than it was. Fungus may be a cause of pulmonary effusion in patients with compromised immunologic defenses.

## Cytology

A cytologic study is performed to detect tumor cells in patients with malignant effusions. Breast and lung are the two most common tumors; lymphoma is the third.

## Carcinoembryonic antigen

Pleural fluid CEA levels are elevated in various malignant (gastrointestinal, breast) conditions. (see p. 194).

## Special tests

The pH of pleural fluid is usually 7.4 or greater. The pH is typically less than 7.2 when empyema is present. The pH may be 7.2 to 7.4 in TB, or malignancy.

In some instances, the rheumatoid factor (see p. 790) and the complement levels (see p. 266) are also measured in pleural fluid.

Pleural fluid antinuclear antibody (ANA) levels and pleural fluid/serum ANA ratios are often used to evaluate pleural effusion secondary to systemic lupus erythematosus.

## Contraindications

- Patients with significant thrombocytopenia

## Potential complications

- Pneumothorax because of puncture of the visceral pleura or entry of air into the pleural space
- Interpleural bleeding because of puncture of tissue or a blood vessel
- Hemoptysis caused by needle puncture of a pulmonary vessel or by inflammation
- Reflex bradycardia and hypotension
- Pulmonary edema
- Seeding of the needle track with tumor when malignant pleural effusion exists

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient.

- Obtain informed consent for this procedure.

**PT** Tell the patient that no fasting or sedation is necessary.

**PT** Inform the patient that movement or coughing should be minimized to avoid inadvertent needle damage to the lung or pleura during the procedure.

- Administer a cough suppressant before the procedure if the patient has a troublesome cough.

### During

- Note the following procedural steps:
    1. The patient is usually placed in an upright position, with the arms and shoulders raised and supported on a padded overhead table. This position spreads the ribs and enlarges the intercostal space for insertion of the needle.
    2. The thoracentesis is performed under strict sterile technique.
    3. The needle insertion site, which is determined by percussion, auscultation, and examination of a chest x-ray image, ultrasound scan, or fluoroscopy, is aseptically cleansed and anesthetized locally.
    4. The needle is positioned in the pleural space, and the fluid is withdrawn with a syringe and a three-way stopcock.
    5. A short polyethylene catheter may be inserted into the pleural space for fluid aspiration; this decreases the risk of puncturing the visceral pleura and inducing a pneumothorax.
  - Monitor the patient's pulse for reflex bradycardia and evaluate the patient for diaphoresis and a feeling of faintness.
  - Note that this procedure is performed by a physician at the patient's bedside, in a procedure room, or in the physician's office in less than 30 minutes.
- PT** Although local anesthetics eliminate pain at the insertion site, tell the patient that he or she may feel a pressure-like pain when the pleura is entered and the fluid is removed.

### After

- Place a small bandage over the needle site. Usually turn the patient on the unaffected side for 1 hour to allow the pleural puncture site to heal.
- Label the specimen with the patient's name, date, source of fluid, and diagnosis. Send the specimen promptly to the laboratory.
- Obtain a chest x-ray study as indicated to check for pneumothorax.
- Monitor the patient's vital signs.
- Observe the patient for coughing or expectoration of blood (hemoptysis), which may indicate trauma to the lung.
- Evaluate the patient for signs and symptoms of pneumothorax, tension pneumothorax, subcutaneous emphysema, and pyogenic infection (e.g., tachypnea, dyspnea, diminished breath sounds, anxiety, restlessness, fever).
- Assess the patient's lung sounds for diminished breath sounds, which could be a sign of pneumothorax.



**Abnormal findings****Exudate**

Empyema  
Pneumonia  
TB, effusion  
Pancreatitis  
Ruptured esophagus  
Tumors  
Lymphoma  
Pulmonary infarction  
Collagen vascular disease  
Drug hypersensitivity

**Transudate**

Cirrhosis  
Congestive heart failure  
Nephrotic syndrome  
Hypoproteinemia  
Trauma

**notes**

## thoracoscopy

**Type of test** Endoscopy

### Normal findings

Normal pleura and lung

### Test explanation and related physiology

This procedure is used to directly visualize the pleura, lung, and mediastinum. Tissue can be obtained for testing. It is also helpful in assisting in the staging and dissection of lung cancers.

With this technique, the parietal pleura, visceral pleura, and mediastinum can be directly visualized. Tumors involving the chest cavity can be staged by direct visualization. Any abnormality can be biopsied. Collections of fluid can be drained and aspirated for testing. Dissection for lung resection can be carried out with the thoracoscope (*video-assisted thoracotomy [VAT]*), thereby minimizing the extent of a thoracotomy incision.

### Contraindications

- Patients with previous lung surgery because it is difficult to obtain access to the free pleural space

### Potential complications

- Bleeding
- Infection or empyema
- Prolonged pneumothorax

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Ensure that informed consent for this procedure is obtained.
- Inform the patient that an open thoracotomy may be required because of the possibility of intrathoracic injury.
- Because the procedure is usually performed with the patient under general anesthesia, follow routine general anesthesia precautions.
- Shave and prepare the patient's chest as ordered.
- Keep the patient NPO (nothing by mouth) after midnight on the day of the test. Intravenous (IV) fluids may be given.

#### During

- Note the following procedural steps:
  1. Thoracoscopy is performed in the operating room. The patient is initially placed in the lateral decubitus position.

2. After the thorax is cleansed, a blunt-tipped (Verres) needle is inserted through a small incision, and the lung is collapsed.
3. A thoracoscope is inserted through a trocar to examine the chest cavity. Other trocars can be placed as conduits for other instrumentation.
4. After the desired procedure is completed, the scope and trocars are removed.
5. Usually a small chest tube is placed to ensure full reexpansion of the lung.
6. The incision(s) is (are) closed with a few skin stitches and covered with an adhesive bandage.

### After

- Assess the patient frequently for signs of bleeding (increased pulse rate, decreased blood pressure). Report any significant findings to the physician.
- Provide analgesics to relieve the minor to moderate pain that may be experienced.
- If a surgical procedure has been performed thoracoscopically, provide appropriate specific postsurgical care.
- A chest x-ray is performed after the procedure to ensure complete reexpansion of the lung.

### Abnormal findings

Primary lung cancer

Metastatic cancer to the lung or pleura

Empyema

Pleural tumor

Pleural infection

Pleural inflammation

Pulmonary infection

### notes

## throat and nose cultures

**Type of test** Microscopic examination

### Normal findings

Negative

### Test explanation and related physiology

Because the throat and nose are normally colonized by many organisms, cultures of these areas serve only to isolate and identify a few particular pathogens (e.g., streptococci, meningococci, gonococci, *Bordetella pertussis*, *Corynebacterium diphtheriae*).

Streptococci are most often sought on a throat culture because beta-hemolytic streptococcal pharyngitis may be followed by rheumatic heart disease or glomerulonephritis. This type of streptococcal infection most frequently affects children between the ages of 3 and 15 years. Throat cultures in adults are indicated only when the patient has a severe or recurrent sore throat often associated with fever and palpable lymphadenopathy.

Rapid immunologic tests (*strep screens*) with antiserum against group A streptococcus antigen are now available and are very accurate. With these kits, the streptococcus organism can be identified directly from the swab specimen without culture. These tests can be performed in about 15 minutes. If the test result is negative, no streptococcus infection exists.

Nasal and nasopharyngeal cultures are often done to screen for infections and carrier states caused by various other organisms, such as *Staphylococcus aureus*, *Haemophilus influenzae*, *Neisseria meningitidis*, respiratory syncytial virus (RSV), and viruses causing rhinitis. Health-care workers in the operating room and newborn nursery may have these cultures to screen potential sources of spread when an outbreak occurs in a hospital setting. These cultures are also used to detect infection in elderly and debilitated patients.

All cultures should be performed before antibiotic therapy is initiated. Otherwise, the antibiotic may interrupt the growth of the organism in the laboratory. Most organisms take approximately 24 hours to grow in the laboratory, and a preliminary report can be given at that time. Occasionally, 48 to 72 hours is required for growth and identification of the organism. Cultures may be repeated on completion of appropriate antibiotic therapy to identify resolution of the infection.

## Interfering factors

- ❖ Drugs that may affect test results include antibiotics and anti-septic mouthwashes.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.

### During

- Obtain a *throat culture* by depressing the tongue with a wooden tongue blade and touching the posterior wall of the throat and areas of inflammation, exudation, or ulceration with a sterile cotton swab. Two swabs are preferred. Growth of streptococcus from both swabs is more accurate, and the second swab can also be used in the strep screen. Avoid touching any other part of the mouth. Place the swabs in a sterile container.
- Obtain a *nasal culture* by gently raising the tip of the nose and inserting a flexible swab into the nares. Rotate the swab against the side of the nares. Remove the swab and place it in an appropriate culture tube.
- Obtain a *nasopharyngeal culture* by gently raising the tip of the nose and inserting a flexible swab along the bottom of the nares. Guide this swab until it reaches the posterior pharynx. Rotate the swab to obtain secretions and then remove it. Place the swab in an appropriate culture tube.
- Wear gloves and handle the specimen as if it were capable of transmitting disease.
- Indicate on the laboratory slip any medications that the patient may be taking that could affect test results.

### After

- Label all specimens and send them immediately to the microbiology laboratory.
- Notify the physician of any positive results so that appropriate antibiotic therapy can be initiated.

## Abnormal findings

Bacterial pathogens (e.g., streptococci)

Respiratory syncytial virus

*H. influenzae* bacteria

## notes


**thromboelastography** (Thromboelastometry)**Type of test** Blood**Normal findings**5.3-12.4 dynes/cm<sup>2</sup>**Possible critical value**> 12.4 dynes/cm<sup>2</sup>**Test explanation and related physiology**

Hemostasis is a balanced process in which the blood forms localized clots when the integrity of the vascular system is breached. Trauma, infection, and inflammation all activate the blood's clotting system, which depends on the interaction of two separate systems: enzymatic proteins (clotting factors, intrinsic and extrinsic systems; see [Figure 41](#)) and platelets. The two systems work in concert to plug defects in the broken vessels. The clots that form in this process need to be of sufficient strength to resist dislodgement. If a particular clotting factor is dysfunctional or absent, as in hemophilia, an insufficient amount of fibrin forms. Similarly, massive consumption of clotting factors in a trauma situation decreases the amount of fibrin formed. Inadequate numbers of platelets resulting from trauma, surgery, or chemotherapy also decrease platelet aggregation, as do genetic disorders, uremia, or medication therapy. Ultimately, reduced fibrin formation or platelet aggregation results in clots of inadequate tensile strength. This *hypocoagulable* state is associated with excessive bleeding. Conversely, endothelial injury, stasis, cancer, genetic diseases, or other *hypercoagulable* states lead to thrombosis formation, causing thromboembolic events.

This test is used to identify patients who are hypercoagulable and may experience a thromboembolic phenomenon when immobile (e.g., after surgery or trauma). It is particularly helpful in cardiac surgery and liver transplantation. It is also used to determine hyperfibrinolysis. This test can evaluate platelet function and determine the percent of platelet inhibition instigated by drugs such as heparin, aspirin, and antiplatelet drugs (Plavix, Ticlid).

**Interfering factors**

- ❖ Drugs that may cause *decreased* thromboelastography include antiplatelet drugs (e.g., ticlopidine), some antibiotics, aspirin, beta-blockers, clofibrate, dextran, ethanol, heparin, nonsteroidal

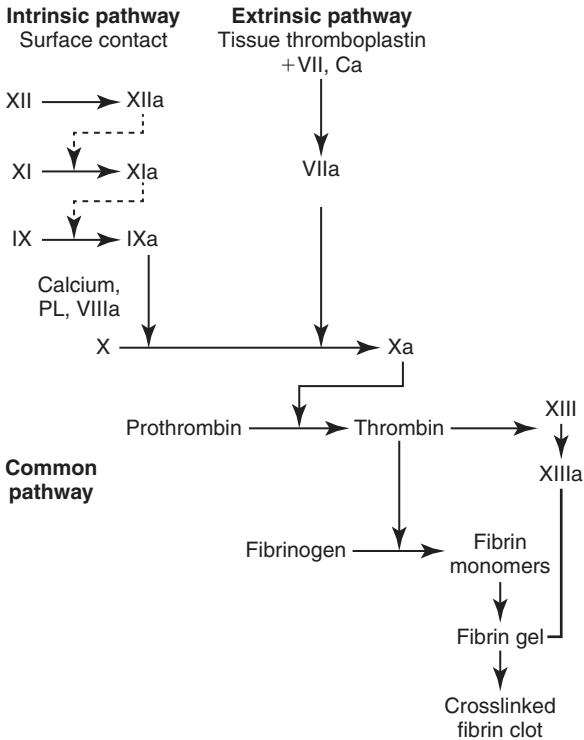

**FIG. 41** Simplified enzymatic cascade of fibrin clot formation.

antiinflammatory drugs, phenothiazines, theophylline, tricyclics, and warfarin sodium (Coumadin).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- If the patient is receiving any drugs that may interfere with normal coagulation or has any diseases such as jaundice, hyperlipidemia, or hemolysis, this should be listed on the laboratory request slip.
- Immediately transfer the specimen to the laboratory.
- Remember that abnormalities in platelet aggregation can prolong bleeding time, and a significant hematoma at the venipuncture site may occur.

## Abnormal findings

### Hypocoagulability

- Factor deficiency
- Anticoagulation
- Thrombocytopenia
- Platelet function abnormalities
- Increased fibrinolysis

### Hypercoagulability

- Factor V Leiden abnormality
- Protein S/C abnormality
- Genetic hypercoagulability
- Idiopathic hypercoagulability

### notes



**thrombosis indicators****fibrin monomers** (Fibrin degradation products [FDPs], Fibrin split products [FSPs])**fibrinopeptide A (FPA)****prothrombin fragment (F1 + 2)****Type of test** Blood**Normal findings**

FDP: &lt; 10 mcg/mL or &lt; 10 mg/L (SI units)

FPA

Male: 0.4-2.6 mg/mL

Female: 0.7-3.1 mg/mL

F1+2: 7.4-103 mcg/L or 0.2-2.8 nmol/L

**Possible critical values**

FDP: &gt; 40 mcg/mL

**Test explanation and related physiology**

Identification of FDPs, FPA, and F1 + 2 is mostly used to document that fibrin clot formation and therefore thrombosis is occurring in patients. These tests support the diagnosis of disseminated intravascular coagulation (DIC). The D-dimer test (see p. 319) is more commonly used to identify DIC or other forms of thrombosis. These tests also provide an indication of the effectiveness of anticoagulation therapy. Finally, they are used to support the diagnosis and follow treatment for hypercoagulable states.

F1 + 2 is liberated when prothrombin is converted to thrombin in reaction 4 of secondary hemostasis (see [Figure 10](#), p. 251). FPA is released into the bloodstream from alpha and beta chains of fibrinogen during its conversion to fibrin.

Measurement of FDPs is an indicator of the activity of the fibrinolytic system. Whenever a fibrin clot or primarily fibrinogen degenerates, fragment monomers called FDPs (X, D, E, and Y) result. These, therefore, are indirect evidence of thrombosis and/or DIC. FDPs also are increased with other secondary fibrinolytic disorders. Thrombolytic therapy used to treat vascular thrombosis is associated with increased FDPs. Streptokinase or urokinase stimulates the conversion of plasminogen to plasmin.

These products of hemostasis and fibrinolysis also may be elevated in patients with extensive malignancy, tissue necrosis, trauma surgery, or gram-negative sepsis. For discussion of D-dimer fibrin degradation products, see p. 319.

## 878 thrombosis indicators

### Interfering factors

- Traumatic venipunctures may increase FPA levels.
- Surgery or massive trauma is associated with increased levels.
- Menstruation may be associated with increased FDP levels.
- ▼ Drugs that may cause *increased* levels include barbiturates, streptokinase, and urokinase.
- ▼ Drugs that may cause *decreased* levels include warfarin and other oral anticoagulants.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: verify with laboratory
- Avoid prolonged use of a tourniquet.
- Draw the sample before initiating heparin therapy.
- Avoid excessive agitation of the blood sample.

### Abnormal findings

#### ▲ Increased levels

DIC  
Heart or vascular surgery  
Thromboembolism  
Thrombosis  
Advanced malignancy  
Severe inflammation  
Postoperative states  
Massive trauma  
Deficiency in proteins S  
and C  
Antithrombin III deficiency

#### ▼ Decreased levels

Anticoagulation therapy

notes

## thyroglobulin (Tg, Thyrogen-stimulated thyroglobulin)

**Type of test** Blood

### Normal findings

Age	Male, ng/mL	Female, ng/mL
0-11 months	0.6-5.5	0.5-5.5
1-11 years	0.6-50.1	0.5-52.1
12 years and older	0.5-53.0	0.5-43.0

### Test explanation and related physiology

Thyroglobulin is the protein precursor of thyroid hormone and is made by both normal, well-differentiated, benign thyroid cells and thyroid cancer cells. Because thyroglobulin is normally only made by thyroid cells, it serves as a useful readout for the presence or absence of thyroid cells, especially after thyroid cancer surgery. This test is primarily used as a tumor marker for well-differentiated thyroid cancers (see p. 194, Tumor markers).

In the treatment of well-differentiated thyroid cancers, it is important to remove as much thyroid tissue as possible so that adjunctive radioactive iodine treatment will not go to residual thyroid gland tissue in the neck but will go instead to any metastatic thyroid cells. If postoperative Tg levels are low, very little thyroid tissue remains.

After thyroidectomy, thyroid hormone replacement is required for normal metabolic function. Thyroid-stimulating hormone (TSH) levels are usually very low when thyroid hormone is replaced adequately. Endogenous stimulation of any residual thyroid cells is minimal in these patients. As a result, Tg and endogenous thyroid hormones are low.

*Thyrogen-stimulated testing* has eliminated the need for withdrawal of thyroid hormone medications and provides a safe and effective method to elevate TSH levels so that even minimal levels of Tg can be detected. This allows patients to undergo periodic thyroid cancer follow-up evaluations while avoiding the often debilitating side effects of hypothyroidism, which is caused by withdrawal of hormone medication.

Thyrogen is a highly purified recombinant source of human TSH. Thyrogen raises serum TSH levels and thereby stimulates Tg production. Normal thyroid-remnant and well-differentiated thyroid tumors display a greater (> 10-fold) serum Tg response to TSH stimulation. If, after thyroid surgery, Thyrogen-stimulated

## 880 thyroglobulin

Tg levels are elevated, either a significant amount of normal thyroid gland was left in the neck or metastatic disease exists.

Thyrogen stimulation is also used for patients undergoing I<sup>131</sup> whole-body scanning for metastatic thyroid cancer. Now, with the use of Thyrogen, the ill effects of hormone withdrawal are not experienced.

### Interfering factors

- Tg levels are decreased in less well-differentiated thyroid cancers.
- Thyrogen stimulation of Tg levels is less in patients whose tumors do not have TSH receptors or whose tumors cannot make Tg.
- Tg autoantibodies cause either underestimation or overestimation of serum Tg measurements made by immunometric assay (IMA) or radioimmunoassay (RIA) methods, respectively.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: serum separator
- Determine whether the patient is to have a whole-body nuclear scan along with the Tg blood test.
- If Thyrogen stimulation is to be used:
  1. Administer Thyrogen intramuscularly to the buttock every 24 hours for two or three doses as ordered.
  2. Collect a venous blood sample in a gold-top (serum separator) tube after 3 days.
- For radioiodine imaging:
  1. The nuclear medicine technologist will administer the radioiodine 24 hours after the final Thyrogen injection.
  2. Scanning is usually performed 48 hours after radioiodine administration. Whole-body images are acquired for a minimum of 30 minutes and/or should contain a minimum of 140,000 counts.
  3. Scanning times for single (spot) images of body regions may be obtained.

### Abnormal findings

#### ▲ Increased levels

Residual thyroid tissue in the neck

Metastatic thyroid cancer

## thyroid fine needle aspiration biopsy (FNAB, Skinny-needle thyroid biopsy, Fine-needle aspiration [FNA])

**Type of test** Microscopic examination of tissue

### Normal findings

Benign

### Test explanation and related physiology

FNAB is used to obtain tissue to rule out cancer in a “cold” (does not light up on a thyroid scan) thyroid nodule or cyst of a patient whose thyroid function is normal. In FNAB, samples of thyroid tissue are obtained by placing a very thin needle into the thyroid nodule and obtaining small pieces of thyroid tissue that are then examined microscopically. The Bethesda System for Reporting Thyroid Cytopathology is used to provide the results of the biopsy. Results are indicated as the following: non-diagnostic (insufficient tissue for evaluation), benign, atypia of undetermined significance, suspicious for a follicular neoplasm, suspicious for malignancy, or malignant. The implied risk of malignancy for each diagnostic category is shown in the following table.

Diagnosis	Chance of malignancy with surgical biopsy, %
Nondiagnostic	1-4
Benign	0-3
Atypia	5-15
Follicular neoplasm	15-30
Suspicious	60-75
Malignant	> 97

FNAB samples can be challenging to interpret and produce indeterminate results in 15% to 30% of cases. Approximately 70% to 80% of the time, the nodules prove to be benign for cancer by surgical pathology.

### Contraindications

- Patients with coagulation disorders because of the risk of excessive bleeding
- Patients with hyperthyroidism, because the needle insertion may instigate thyroid storm and toxic nodular goiters are not malignant.

## 882 thyroid fine needle aspiration biopsy

### Potential complications

- Hemorrhage from the highly vascular thyroid tissue
- Cyst formation in the thyroid gland
- Infection when an open biopsy is performed

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

- Ensure that the physician obtains informed consent for this procedure.

**PT** Inform the patient that no abstinence from food or drink is required.

**PT** Tell the patient that no sedative is required.

- Note that the needle stick may be done at the bedside.
- If ultrasound guidance is to be used, note that the needle stick is performed in the radiology or ultrasonography department.

#### During

- Note the following procedural steps:
  1. The patient is placed in a supine position with a sandbag or pillow under the shoulder to hyperextend the neck.
  2. Under sterile conditions, the skin overlying the thyroid is infiltrated with a local anesthetic (lidocaine).
  3. If the nodule can be felt, a biopsy can be performed in the doctor's office. When the nodule is not palpable, ultrasound is used to help guide the biopsy.
  4. The patient holds her or his breath while the needle is rocked gently to obtain as much tissue as possible.
  5. The needle is then withdrawn and tissue is placed on a glass slide.

#### After

- Apply pressure over the thyroid area to minimize bleeding.
- Note that a physician performs this procedure in approximately 10 minutes.

### Abnormal findings

Malignancy

notes

## thyroid scanning (Thyroid scintiscan)

**Type of test** Nuclear scan

### Normal findings

Normal size, shape, position, and function of the thyroid gland  
No areas of decreased or increased uptake

### Test explanation and related physiology

Thyroid scanning allows the size, shape, position, and physiologic function of the thyroid gland to be determined with the use of radionuclear scanning.

Thyroid nodules are easily detected by this technique. Nodules are classified as functioning (warm/hot) or nonfunctioning (cold), depending on the amount of radionuclide taken up by the nodule. A functioning nodule could represent a benign adenoma or a localized toxic goiter. A nonfunctioning nodule may represent a cyst, carcinoma, nonfunctioning adenoma or goiter, lymphoma, or localized area of thyroiditis.

Scanning is useful in patients with

- A neck or substernal mass
- A thyroid nodule. Thyroid cancers are usually nonfunctioning (cold) nodules.
- Hyperthyroidism. Scanning assists in differentiating Graves disease (diffusely enlarged hyperfunctioning thyroid gland) from Plummer disease (nodular hyperfunctioning gland).
- Metastatic tumors without a known primary site. A normal scan excludes the thyroid gland as a possible primary site.
- A well-differentiated form of thyroid cancer. Areas of metastasis may show up on subsequent whole-body nuclear scans.

Thyroid scanning is usually preceded by a thyroid uptake scan. An iodine-123 capsule is given to the patient 6 to 24 hours before measuring iodine uptake. After uptake counts are obtained, Tc 99m<sup>04</sup> is administered intravenously. Scan images are then obtained as described earlier.

Another form of thyroid scan is called the *whole-body thyroid scan*. This scan is performed on patients who have previously had a thyroid cancer treated. An iodine-131 capsule or solution, Tc 99m<sup>04</sup>, or iodine-123 is given orally, and the entire body is scanned to look for metastatic thyroid tissue. A hot spot indicates recurrent tumor. This test is routinely performed (every 1-2 years) on patients who have had a thyroid cancer larger than 1 cm.

## Contraindications

- Patients who are allergic to iodine
- Patients who are pregnant unless the benefits outweigh the risks

## Potential complications

- Radiation-induced oncogenesis. This complication is minimized if technetium or low-radioactive iodine isomers are used instead of iodine-131.

## Interfering factors

- Iodine-containing foods
- Recent administration of x-ray contrast agents
- ✚ Drugs that may affect test results include cough medicines, multiple vitamins, oral contraceptives, and thyroid drugs.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- Check the patient for allergies to iodine.
- PT** Instruct the patient about medications that need to be restricted for weeks before the test (e.g., thyroid drugs, medications containing iodine).
- Obtain a history concerning previous contrast x-ray studies, nuclear scanning, or intake of any thyroid-suppressive or anti-thyroid drugs.
- PT** Tell the patient that fasting is usually not required. Check with the laboratory.

### During

- Note the following procedural steps:
  1. A standard dose of iodine-123 is usually given to the patient by mouth 6 to 24 hours before scanning. The capsule or solution is tasteless.
  2. Iodine uptake is measured. IV technetium is administered, and thyroid imaging is performed 2 hours later.
  3. At the designated time, the patient is placed in a supine position, and anterolateral images of the thyroid area are obtained. The whole-body imaging protocol is somewhat different.
  4. The radioactive counts are recorded and displayed.



- Note that this study is performed by a radiologic technologist in less than 30 minutes and is then interpreted by a nuclear medicine physician.
- PT** Tell the patient that no discomfort is associated with this study.

### After

- PT** Usually the dose of radioactivity used in this test is minimal and considered harmless. No isolation and no special urine precautions are needed. However, if higher doses of radionuclide are used, isolation for 24 hours may be recommended.

### Abnormal findings

Adenoma

Toxic and nontoxic goiter

Cyst

Carcinoma

Lymphoma

Thyroiditis

Graves disease

Plummer disease

Metastasis

Hyperthyroidism

Hypothyroidism

Hashimoto thyroiditis

**notes**

**thyroid-stimulating hormone** (TSH, Thyrotropin, TRH stimulation test)**Type of test** Blood**Normal findings**Adult: 2-10  $\mu\text{U}/\text{mL}$  or 2-10  $\text{mU}/\text{L}$  (SI units)Newborn: 3-18  $\mu\text{U}/\text{mL}$  or 3-18  $\text{mU}/\text{L}$ Cord: 3-12  $\mu\text{U}/\text{mL}$  or 3-12  $\text{mU}/\text{L}$ 

(Values vary among laboratories.)

**Test explanation and related physiology**

The TSH concentration aids in differentiating primary from secondary hypothyroidism. Pituitary TSH secretion is stimulated by hypothalamic thyroid-releasing hormone (TRH). Low levels of triiodothyronine and thyroxine ( $\text{T}_3$  and  $\text{T}_4$ ) are the underlying stimuli for TRH and TSH. Therefore a compensatory elevation of TRH and TSH occurs in patients with primary hypothyroid states such as surgical or radioactive thyroid ablation; patients with burned-out thyroiditis, thyroid agenesis, idiopathic hypothyroidism, or congenital hypothyroidism; and patients taking antithyroid medications.

In secondary hypothyroidism, the function of the hypothalamus or pituitary gland is faulty because of tumor, trauma, or infarction. Thus TRH and TSH cannot be secreted, and plasma levels of these hormones are near 0 despite low  $\text{T}_3$  and  $\text{T}_4$  levels.

The *TRH stimulation test* is sometimes used to stimulate low levels of TSH to identify primary from secondary hypothyroidism in cases in which TSH is low. However, this test is not commonly used because extremely low levels of TSH can be identified now with the use of immunoassays.

The TSH test is used as well to monitor exogenous thyroid replacement. The goal of thyroid replacement therapy is to provide an adequate amount of thyroid medication so that TSH secretion is in the low normal range, indicating a euthyroid state. This test is also done to detect primary hypothyroidism in newborns with low screening  $\text{T}_4$  levels. TSH and  $\text{T}_4$  levels are frequently measured to differentiate pituitary from thyroid dysfunction. A decreased  $\text{T}_4$  with a normal or elevated TSH level can indicate a thyroid disorder. A decreased  $\text{T}_4$  with a decreased TSH level can indicate a pituitary disorder.

**Interfering factors**

- Recent radioisotope administration may affect test results.

- Severe illness may cause decreased TSH levels.
- ✚ Drugs that may cause *increased* levels include antithyroid medications, lithium, potassium iodide, and TSH injection.
- ✚ Drugs that may cause *decreased* levels include aspirin, dopamine, heparin, steroids, and T<sub>3</sub>.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Use a heel stick to obtain blood from newborns.

### Abnormal findings

#### ▲ Increased levels

Primary hypothyroidism  
 (thyroid dysfunction)  
 Thyroiditis  
 Thyroid agenesis  
 Congenital hypothyroidism  
 Large doses of iodine  
 Severe and chronic illnesses  
 Pituitary TSH-secreting  
 tumor

#### ▼ Decreased levels

Secondary  
 hypothyroidism  
 (pituitary dysfunction)  
 Hyperthyroidism  
 Pituitary hypofunction

notes

**thyroid-stimulating hormone stimulation test (TSH stimulation test)****Type of test** Blood**Normal findings**

Increased thyroid function with administration of exogenous TSH

**Test explanation and related physiology**

The TSH stimulation test is used to differentiate *primary* (or thyroidal) hypothyroidism from *secondary* (or hypothalamic-pituitary) hypothyroidism. Normal people and patients with hypothalamic-pituitary hypothyroidism can increase thyroid function when exogenous TSH is given. However, patients with primary thyroidal hypothyroidism do not; their thyroid glands are inadequate and cannot function no matter how much stimulation they receive. Patients with less than a 10% increase in radioactive iodine uptake (RAIU) or less than a 1.5 mcg/dL rise in  $T_4$  are considered to have a primary cause for their hypothyroid state. If the initially low uptake is caused by inadequate pituitary stimulation of an intrinsically normal thyroid gland, the RAIU should increase at least 10% and the  $T_4$  level should rise 1.5 mcg/dL or more. This is characteristic of secondary hypothyroidism.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: Verify with laboratory.
- Obtain baseline levels of RAIU or  $T_4$  (see p. 894) as indicated.
- Administer the prescribed dose of TSH intramuscularly for 3 days.
- Repeat the levels of RAIU or  $T_4$  as indicated.

**Abnormal findings**

Primary (thyroidal) hypothyroidism

Secondary (hypothalamic-pituitary) hypothyroidism

**notes**

**thyroid-stimulating immunoglobulins** (TSIs, Long-acting thyroid stimulator [LATS], Thyroid-binding inhibitory immunoglobulin [TBII], Thyrotropin receptor antibody, Thyroid-stimulating hormone receptor [TSHR] antibodies)

**Type of test** Blood

### Normal findings

TSI: < 130% of basal activity

TBII: < 10% inhibition

### Test explanation and related physiology

TSIs represent a group of immunoglobulin G (IgG) antibodies directed against the thyroid cell receptor for TSH and are associated with autoimmune thyroid disease states, such as chronic thyroiditis and Graves disease. (These autoantibodies bind and transactivate the TSHRs). This instigates stimulation of the thyroid gland independent of the normal feedback-regulated TSH stimulation. This, in turn, will stimulate the release of thyroid hormones from the thyroid cells. Some patients with Graves disease also have TSHR-blocking antibodies, which do not transactivate the TSHR. The balance between TSI and TSHR-blocking antibodies, as well as their individual titers, is thought to be a determinant of Graves disease severity.

The use of these antibodies is helpful in the evaluation of patients for whom the diagnosis of Graves disease is confused by conflicting data (e.g., subclinical Graves hyperthyroidism or euthyroid patients with ophthalmopathy). In these cases, the antibodies help determine and support the diagnosis of Graves disease.

Because TSIs can cross the placenta, they may be found in neonates whose mothers have Graves disease.

### Interfering factors

- Recent administration of radioactive iodine may affect test results.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or gold
- Notify the laboratory if the patient has received radioactive iodine in the preceding 2 days.
- Handle the blood sample gently. Hemolysis may interfere with interpretation of test results.

## Abnormal findings

### ▲ Increased levels

- Hyperthyroidism
- Malignant exophthalmos
- Graves disease
- Hashimoto thyroiditis
- Neonatal thyrotoxicosis

notes

## thyroid ultrasound (Thyroid echogram, Thyroid sonogram)

**Type of test** Ultrasound

### Normal findings

Normal size, shape, and position of the thyroid gland

### Test explanation and related physiology

Ultrasound examination of the thyroid gland is valuable for distinguishing cystic from solid thyroid nodules. If the nodule is found to be purely cystic (fluid filled), the fluid can simply be aspirated and surgery avoided. If the nodule has a mixed or solid appearance, however, a tumor may be present, and FNA (see p. 881) surgery may be required. This study may be repeated at intervals to determine the response of a thyroid mass to medical therapy.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that breathing or swallowing will not be affected by the placement of a transducer on the neck.
- PT** Inform the patient that a lubricant will be applied to the neck to ensure effective transmission of sound waves.
- PT** Tell the patient that no fasting or sedation is required.

#### During

- Note the following procedural steps:
  1. The patient is taken to the ultrasonography department and placed in the supine position with the neck hyperextended.
  2. Gel is applied to the patient's neck.
  3. A sound transducer is passed over the nodule.
  4. Photographs are taken of the image displayed.
- A technologist usually performs this study in about 15 minutes and a radiologist evaluates the results.
- PT** Tell the patient that no discomfort is associated with this study.

#### After

- Assist the patient in removing the lubricant from the neck.

### Abnormal findings

Cyst  
Tumor  
Thyroid adenoma

Thyroid carcinoma  
Goiter

**thyrotropin-releasing hormone stimulation test** (TRH stimulation test, Thyrotropin-releasing factor stimulation test [TRF stimulation test])**Type of test** Blood**Normal findings**Baseline TSH:  $< 10 \mu\text{U}/\text{mL}$ 

Stimulated TSH: more than double the baseline

**Test explanation and related physiology**

The TRH stimulation test assesses the anterior pituitary gland via its secretion of TSH in response to an IV injection of TRH. After the TRH injection, the normally functioning pituitary gland should secrete TSH. In hyperthyroidism, either a slight or no increase in the TSH level is seen because pituitary TSH production is suppressed by the direct effect of excess circulating  $T_4$  and  $T_3$  on the pituitary gland. A normal result is considered reliable evidence for excluding the diagnosis of thyrotoxicosis. Since the development of very sensitive radioimmunoassay for TSH, the TRH stimulation test is no longer required to diagnose hyperthyroidism. However, it still has a role in the evaluation of pituitary deficiency.

In addition to assessing the responsiveness of the anterior pituitary gland, this test aids in the detection of primary, secondary, and tertiary hypothyroidism. In primary hypothyroidism (thyroid gland failure), the increase in the TSH level is two or more times the normal result. With secondary hypothyroidism (anterior pituitary failure), no TSH response occurs. Tertiary hypothyroidism (hypothalamic failure) may be diagnosed by a delayed rise in the TSH level. Multiple injections of TRH may be needed to induce the appropriate TSH response in this case.

The TRH test also may be useful in differentiating primary depression from manic-depressive psychiatric illness and from secondary types of depression. In primary depression, whereas the TSH response is blunted in most patients, patients with other types of depression have a normal TRH-induced TSH response.

**Interfering factors**

- Pregnancy may increase the TSH response to TRH.
- ✦ Drugs that may modify the TSH response include antithyroid drugs, aspirin, corticosteroids, estrogens, levodopa, and  $T_4$ .



## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: verify with laboratory
- PT** Instruct the patient to discontinue thyroid preparations for 3 to 4 weeks before the TRH test, if indicated.
- Assess the patient for medications currently being taken.
- Administer a prescribed IV bolus of TRH.
- Obtain venous blood samples at intervals, and measure for TSH levels.
- Indicate on the laboratory slip if the patient is pregnant.

## Abnormal findings

Hyperthyroidism

Hypothyroidism

Psychiatric primary depression

Acute starvation

Old age (especially in men)

Pregnancy

## notes

## thyroxine, total and free ( $T_4$ , Thyroxine screen, $FT_4$ )

**Type of test** Blood

### Normal findings

Free  $T_4$ :

0-4 days: 2-6 ng/dL or 26-68 pmol/L (SI units)

2 weeks-20 years: 0.8-2 ng/dL or 10-26 pmol/L (SI units)

Adult: 0.8-2.8 ng/dL or 10-36 pmol/L (SI units)

Total  $T_4$ :

1-3 days: 11-22 mcg/dL or 152-292 nmol/L (SI units)

1-2 weeks: 10-16 mcg/dL or 126-214 nmol/L (SI units)

1-12 months: 8-16 mcg/dL or 101-213 nmol/L (SI units)

1-5 years: 7-15 mcg/dL or 94-194 nmol/L (SI units)

5-10 years: 6-13 mcg/dL or 83-172 nmol/L (SI units)

10-15 years: 5-12 mcg/dL or 72-151 nmol/L (SI units)

Adult male: 4-12 mcg/dL or 59-135 nmol/L (SI units)

Adult female: 5-12 mcg/dL or 71-142 nmol/L (SI units)

Adult > 60 years: 5-11 mcg/dL or 64-142 nmol/L (SI units)

Pregnancy: 9-14 mcg/dL or 117-181 nmol/L (SI units)

### Test explanation and related physiology

The two main hormones secreted by the thyroid gland are thyroxine, which contains four atoms of  $T_4$  and  $T_3$  (see p. 910). More than 90% of thyroid hormone is made up of  $T_4$ . Thyroid hormones circulate primarily bound to carrier proteins (e.g., TBG, albumin, or transthyretin); a small fraction circulates unbound (free). Only the free forms are metabolically active. Ninety-nine percent of  $T_4$  is bound to proteins (TBG, albumin, or transthyretin). Total  $T_4$  measurements consist of both the bound and unbound fractions. Free  $T_4$  is a measure of unbound metabolically active  $T_4$ . Thyroxine tests are used to determine thyroid function. Greater than normal levels indicate hyperthyroid states, and subnormal values are seen in hypothyroid states.  $T_4$  and TSH are used to monitor thyroid replacement and suppressive therapy.

Abnormalities in protein levels can have a significant effect on the results of the total  $T_4$ . Pregnancy and hormone replacement therapy increase TBG and may cause total  $T_4$  to be falsely elevated, suggesting that hyperthyroidism exists when in fact the patient is euthyroid. If the free  $T_4$  is measured in these patients, it would be normal, indicating that free  $T_4$  is a more accurate indicator of thyroid function than total  $T_4$ . In cases in which TBG is reduced (e.g., hypoproteinemia), the total  $T_4$  is likewise reduced,

suggesting hypothyroidism. Measurement of free  $T_4$  would indicate normal levels and thereby discount the abnormal total  $T_4$  as merely a result of the reduced TBG and not as a result of hypothyroidism.

### Interfering factors

- Neonates have higher free  $T_4$  levels than older children and adults.
- Prior use of iodinated radioisotopes or iodinated contrast material can alter test results.
- Pregnancy causes increased total  $T_4$  levels.
- ✚ Drugs that *increase* free  $T_4$  levels include aspirin, danazol, heparin, and propranolol.
- ✚ Drugs that *decrease* free  $T_4$  levels include furosemide, methadone, phenytoins, and rifampicin.
- ✚ Exogenously administered thyroxine causes *increased* free  $T_4$ .
- ✚ Drugs that may cause *increased* total  $T_4$  levels include clofibrate, estrogens, heroin, methadone, and oral contraceptives.
- ✚ Drugs that may cause *decreased*  $T_4$  levels include anabolic steroids, androgens, antithyroid drugs (e.g., propylthiouracil), lithium, phenytoin, and propranolol.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Evaluate the patient's medication history.
- PT** If indicated, instruct the patient to stop exogenous  $T_4$  medication 1 month before testing.
- List on the laboratory slip any drugs that may affect test results.

## Abnormal findings

### ▲ Increased levels

- Graves disease
- Toxic thyroid adenoma
- Acute thyroiditis
- Familial dysalbuminemic hyperthyroxinemia
- Factitious hyperthyroidism
- Struma ovarii
- Pregnancy
- Hepatitis
- Congenital hyperproteinemia
- Thyroid cancer
- Toxic multinodular goiter

### ▼ Decreased levels

- Congenital hypothyroidism
- Surgical ablation of the thyroid
- Myxedema
- Pituitary insufficiency
- Hypothalamic failure
- Iodine insufficiency
- Renal failure
- Cushing syndrome
- Cirrhosis
- Hashimoto thyroiditis
- Thyroid agenesis
- TSH receptor defects

## notes

**thyroxine-binding globulin (TBG, Thyroid-binding globulin)****Type of test** Blood**Normal findings**

Age	Male (mg/dL)	Female (mg/dL)
1-5 days	2.2-5.9	2.2-5.9
1-11 months	3.1-5.6	3-5.6
1-9 years	2.5-5	2.5-5
10-19 years	2.1-4.6	2.1-4.6
> 20 years	1.2-2.5	1.4-3

Oral contraceptive use: 1.5-5.5 mg/dL

Pregnancy (third trimester): 4.7-5.9 mg/dL

**Test explanation and related physiology**

TBG is the major thyroid hormone protein carrier. When it is elevated, total  $T_3$  and total  $T_4$  are also elevated. This may give the false indication that the patient has hyperthyroidism when in fact the patient just has elevated TBGs. When total  $T_4$  is elevated, one must ascertain whether that elevation is caused by an elevation in TBG or a real elevation in  $T_4$  alone, which is associated with hyperthyroidism.

The most common causes of elevated TBGs are pregnancy, hormone replacement therapy, and the use of oral contraceptives. Elevated TBGs are also present in some cases of porphyria and in infectious hepatitis. Decreased TBGs are commonly associated with other causes of hypoproteinemia (e.g., nephrotic syndrome, gastrointestinal malabsorption, and malnutrition).

**Interfering factors**

- ✦ Drugs that *increase* TBG include estrogens, methadone, oral contraceptives, and tamoxifen.
- ✦ Drugs that *decrease* TBG include anabolic steroids, androgens, danazol, phenytoin, propranolol, and steroids.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- List on the laboratory slip any drugs that may affect test results.

### Abnormal findings

#### ▲ Increased levels

- Pregnancy
- Estrogen replacement therapy
- Estrogen-producing tumors
- Infectious hepatitis
- Genetic increased TBG
- Acute intermittent porphyria

#### ▼ Decreased levels

- Protein-losing enteropathy
- Protein-losing nephropathy
- Malnutrition
- Testosterone-producing tumors
- Ovarian failure
- Major stress

#### notes

## TORCH test

The term *TORCH* (toxoplasmosis, other, rubella, cytomegalovirus, herpes) has been applied to maternal infections with recognized detrimental effects on the fetus. TORCH testing refers to the testing for IgG (indicating past infection) and IgM (indicating recent infection) antibodies to the particular infectious agents described. Included in the category of *other* are such infections as syphilis. All of these tests are discussed separately:

- Toxoplasmosis, p. 902
- Rubella, p. 793
- Cytomegalovirus, p. 317
- Herpes virus, p. 498

The direct effects of these infections may be immediate or delayed. A majority of infants infected by toxoplasmosis in utero are asymptomatic at birth with neurologic sequelae appearing later in life. Prenatal rubella infections, on the other hand, may severely damage the fetus, causing congenital heart disease and mental retardation. Congenital cytomegalovirus infections can result in infant death, hearing loss, cerebral palsy, mental retardation, or chorioretinitis. Neonatal herpes infection (can occur in utero or during birth) presents as localized infection to the skin, eye, and mouth (SEM) or as disseminated infection (involving multiple organs, such as the liver, lung, adrenal glands, or brain). The earlier these infections are recognized, the earlier they can be treated or steps can be made to preclude the long-term effects of the disease. Other effects of TORCH infections on the fetus may be indirect and precipitate abortion or premature labor.

The accuracy of prenatal serology testing is variable. TORCH infections can be more accurately diagnosed by direct identification of the organism by polymerase chain reaction (PCR) testing or culture of amniotic fluid.

### notes

**total blood volume** (TBV, Red blood cell [RBC] volume)**Type of test** Nuclear scan**Normal findings**

Approximately 70 mL per kg of body weight

**Test explanation and related physiology**

Measurement of TBV is performed using radionuclides. This measurement is an accurate indicator of true plasma (liquid components of blood) measurement. Based on the patient's height, weight, gender, and body composition, TBV can determine whether the measured volumes are high or low compared with what would be ideal for the particular patient. The report indicates actual volumes for TBV and RBCs that deviate from normal.

Radioiodine-labeled albumin is injected intravenously. Blood is withdrawn every 5 minutes for five samples. The radioactivity is counted and compared with what would be considered normal. A lower amount of radioactivity in the sample indicates a higher plasma volume. The hematocrit is then used to derive the RBC volume. The TBV is obtained by adding the plasma volume and the RBC volume.

This information may be useful in the following clinical circumstances:

- *Congestive heart failure.* The actual amount of fluid overload can be calculated, and diuresis can be more appropriately determined.
- *Presurgery.* The patient's fluid status can be accurately determined as can RBC status.
- *Acutely ill patients.* There are often large fluid shifts in these patients, and TBV may help in guiding IV fluid replacement.
- *Azotemia.* Measurement of TBV will indicate if azotemia is prerenal (hypovolemia) or primary renal.
- *Hypertension.* TBV may indicate plasma volume overload versus vascular constriction.
- *Anemia.* TBV and RBC volumes can indicate accurately the extent of anemia that otherwise could be affected by fluid status, etc.



## Procedure and patient care

### Before

- PT** Explain the procedure and tell the patient that no fasting is required.

### During

- Obtain venous access.
- Fifteen minutes after radionuclide injection, the first venous blood sample is collected in a red-top tube.
- Similar venous blood samples are obtained every 5 minutes for a total of five samples.

### After

- Apply pressure to the IV site upon removal after extracting the last sample.

## Abnormal findings

### ▲ Increased levels

Hypervolemia  
Hypertension  
Congestive heart failure  
Primary renal disease  
Polycythemia vera

### ▼ Decreased levels

Dehydration  
Hypovolemia  
Acute bleeding  
Anemia

notes

## toxoplasmosis antibody titer

**Type of test** Blood

### Normal findings

Ig G titers < 1:16 indicate no previous infection.

Ig G titers 1:16-1:256 are usually prevalent in the general population.

Ig G titers > 1:256 suggest recent infection.

Ig M titers > 1:256 indicate acute infection.

### Test explanation and related physiology

Toxoplasmosis is a protozoan disease caused by *Toxoplasma gondii*, which is found in poorly cooked or raw meat and in cat feces. The Centers for Disease Control and Prevention recommends that patients who are pregnant be serologically tested for this disease.

The presence of antibodies before pregnancy indicates prior exposure and chronic asymptomatic infection. The presence of these antibodies probably ensures protection against congenital toxoplasmosis in the child. Fetal infection occurs if the mother acquires toxoplasmosis after conception and passes it to the fetus through the placenta. Repeat testing of pregnant patients with low or negative titers may be done before the 20th week and before delivery to identify antibody converters and determine appropriate therapy (e.g., therapeutic abortion at 20 weeks, treatment during the remainder of the pregnancy, or treatment of the newborn). Hydrocephaly, microcephaly, chronic retinitis, and seizures are complications of congenital toxoplasmosis.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red

### Abnormal findings

Toxoplasmosis infection

### notes

## transesophageal echocardiography (TEE)

**Type of test** Endoscopy and ultrasonography

### Normal findings

Normal position, size, and movement of the heart muscle, valves, and heart chambers

### Test explanation and related physiology

TEE provides ultrasonic imaging of the heart from a retrocardiac vantage point, avoiding interference with the ultrasound by the interposed subcutaneous tissue, bony thorax, and lungs. In this procedure, a high-frequency ultrasound transducer placed in the esophagus by endoscopy provides better resolution than that of images obtained with routine transthoracic echocardiography (see p. 339). For TEE, the transducer is positioned behind the heart (Figure 42).

TEE is helpful in the evaluation of structures that are inaccessible or poorly visualized by the transthoracic probe approach. It is especially helpful in patients who are obese or have large lung-air spaces.


This test is performed for the following reasons:

- To better visualize the heart valves
- To differentiate intracardiac from extracardiac masses and tumors
- To better visualize the atrial septum (for atrial septal defects)
- To diagnose thoracic aortic dissection
- To better detect valvular vegetation indicative of endocarditis
- To determine cardiac sources of arterial embolism
- To detect coronary artery disease and ischemia

Ischemic muscle movement is much different from normal muscle movement; therefore TEE is a very sensitive indicator of myocardial ischemia. TEE can be used to monitor patients undergoing major abdominal, peripheral vascular, and carotid artery procedures who are at high risk for intraoperative ischemia because of coronary artery disease.

### Contraindications

- Patients with known upper esophageal pathology
- Patients with known esophageal varices
- Patients with Zenker diverticulum
- Patients with esophageal abnormalities (e.g., stricture)
- Patients with bleeding disorders
- Patients who have had prior esophageal surgery
- Patients who cannot cooperate with the procedure

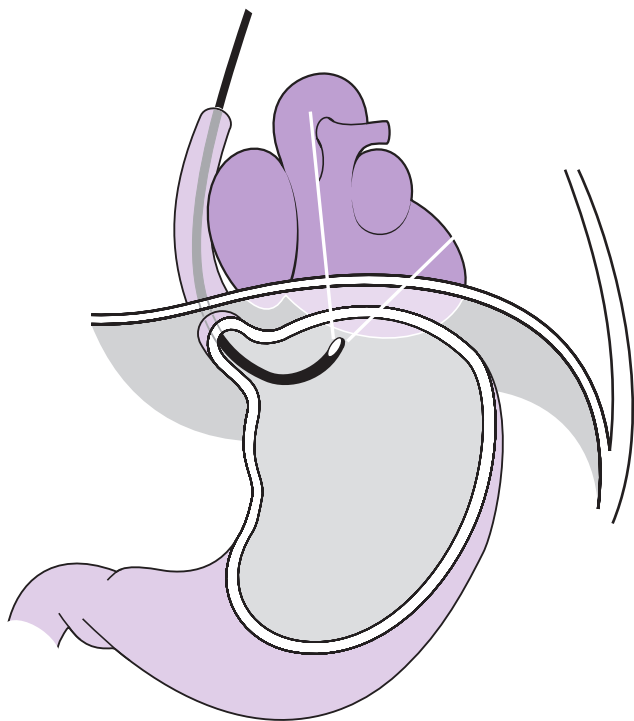

**FIG. 42** Transesophageal echocardiography. The diagram illustrates the location of the transesophageal endoscope in the esophagus.

### Potential complications

- Esophageal perforation or bleeding
- Cardiac arrhythmias

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to fast for 4 to 6 hours before the test.
- PT** Tell the patient to remove all oral prostheses.
- Obtain IV access.

## During

- IV sedation is commonly provided.
- Apply electrocardiography leads and continually monitor heart rhythm.
- Apply a blood pressure cuff, and monitor blood pressure.
- Pulse oximetry is monitored to determine oxygen saturation.
- Note the following procedural steps:
  1. The pharynx is anesthetized with a local topical agent to depress the gag reflex.
  2. The patient is placed in the left lateral decubitus position.
  3. The transducer is advanced into position behind the heart.
  4. Ultrasound images are displayed on a monitor.
- The procedure is performed by a cardiologist and/or a gastrointestinal endoscopist in approximately 20 minutes in the endoscopy suite. It also can be performed at the bedside.
- Very little discomfort is associated with this test.

## After

- Observe the patient closely for approximately 1 hour after the procedure until the effects of sedation have worn off.

## Abnormal findings

Myocardial ischemia

Myocardial infarction

Valvular heart disease

Intracardiac thrombi

Cardiac valvular vegetation

Ventricular and atrial septal defects

Cardiomyopathy

Marked cardiac chamber dilation

Cardiac tumors

Aortic aneurysm or dissection

Aortic plaque

Pulmonary hypertension

Anomalous pulmonary veins

## notes

**transferrin receptor (TfR) assay****Type of test** Blood**Normal findings**

Men: 2-5 mg/L

Women: 1.9-4.4 mg/L

(Results vary depending on the testing method.)

**Test explanation and related physiology**

Both iron metabolism and transport are altered in chronic illness. Differentiation of the anemia of chronic disease (anemia of inflammation or anemia of aging) from iron deficiency anemia is by measurement of the serum TfR concentration.

TfR assists with absorption of iron into cells. TfR is increased when erythropoiesis is enhanced (such as often occurs in iron deficiency). The concentration of cell surface TfR is carefully regulated by TfR mRNA, according to the internal iron content of the cell and its individual iron requirements. Whereas iron-deficient cells contain increased numbers of receptors, receptor numbers are downregulated in iron-replete cells.

An increased mean TfR concentration is noted in patients with iron-deficiency anemia compared with patients with anemia secondary to chronic disease. Calculation of the ratio of TfR/log ferritin concentration provides an even higher sensitivity and specificity for the detection of Fe deficiency.

TfR is also useful in distinguishing iron-deficiency anemia from situations that are commonly encountered in childhood, adolescence, and during pregnancy when iron stores are uniformly low to absent. In these situations, iron-deficient erythropoiesis is not necessarily present, and TfR levels are not elevated. Finally, in situations in which iron-deficiency anemia coexists with anemia of chronic disease, TfR concentrations increase secondary to the underlying iron deficiency, thus avoiding the need for a bone marrow examination.

**Interfering factors**

- Individuals who live at high altitudes have a reference range that extends 6% higher than the upper level of this reference interval.
- Results are related to ethnicity. Individuals of African descent can be expected to have higher levels.
- ✚ Drugs that may cause *increased* TfR levels include recombinant human erythropoietins.

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or green

## Abnormal findings

### ▲ Increased plasma TfR levels

Iron-deficiency anemia  
Reticulocytosis  
Erythropoietin therapy

### ▼ Decreased plasma TfR levels

Hemochromatosis

notes

## triglycerides (TGs)

**Type of test** Blood**Normal findings**

Adult/elderly

Male: 40-160 mg/dL or 0.45-1.81 mmol/L (SI units)

Female: 35-135 mg/dL or 0.40-1.52 mmol/L (SI units)

Child/adolescent

Age, years	Male, mg/dL	Female, mg/dL
0-5	30-86	32-99
6-11	31-108	35-114
12-15	36-138	41-138
16-19	40-163	40-128

**Possible critical values**

&gt; 400 mg/dL

**Test explanation and related physiology**

TGs are a form of fat that exists in the bloodstream. They are transported by very-low-density lipoproteins (VLDLs) and low-density lipoproteins (LDLs). TGs are produced in the liver by using glycerol and other fatty acids as building blocks. TGs act as a storage source for energy. When TG levels in the blood are in excess, TGs are deposited into the fatty tissues. TGs are a part of a lipid profile that also evaluates cholesterol (see p. 235) and lipoproteins (see p. 565). A lipid profile is performed to assess the risk of coronary and vascular disease.

**Interfering factors**

- Ingestion of fatty meals may cause increased TG levels.
- Ingestion of alcohol may cause increased levels.
- Pregnancy may cause increased levels.
- ✚ Drugs that may cause *increased* TG levels include cholestyramine, estrogens, and oral contraceptives.
- ✚ Drugs that may cause *decreased* levels include ascorbic acid, asparaginase, clofibrate, colestipol, fibrates, and statins.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: yes (12-14 hours)
- Blood tube commonly used: red



- PT** Tell the patient not to drink alcohol for 24 hours before the test.
- PT** Inform the patient that dietary indiscretion for as long as 2 weeks before this test will influence results.
  - Mark the patient's age and gender on the laboratory slip.
- PT** Instruct the patient with increased TG levels regarding diet, exercise, and appropriate weight.

### Abnormal findings

#### ▲ Increased levels

Glycogen storage disease  
Hyperlipidemias  
Hypothyroidism  
High-carbohydrate diet  
Poorly controlled diabetes  
Risk of arteriosclerotic  
occlusive coronary disease  
and peripheral vascular  
disease  
Nephrotic syndrome  
Hypertension  
Alcoholic cirrhosis  
Pregnancy  
Myocardial infarction

#### ▼ Decreased levels

Malabsorption syndrome  
Malnutrition  
Hyperthyroidism

notes

**triiodothyronine** ( $T_3$  radioimmunoassay)**Type of test** Blood**Normal findings**

1-3 days: 100-740 ng/dL

1-11 months: 105-245 ng/dL

1-5 years: 105-270 ng/dL

6-10 years: 95-240 ng/dL

11-15 years: 80-215 ng/dL

16-20 years: 80-210 ng/dL

20-50 years: 70-205 ng/dL or 1.2-3.4 nmol/L (SI units)

&gt;50 years: 40-180 ng/dL or 0.6-2.8 nmol/L (SI units)

Reversed  $T_3$ : 10-24 ng/dL**Test explanation and related physiology**

As with the thyroxine ( $T_4$ ) test, the serum  $T_3$  test is an accurate measure of thyroid function.  $T_3$  is less stable than  $T_4$  and occurs in minute quantities in the active form. Only about 7% to 10% of thyroid hormone is composed of  $T_3$ . Seventy percent of that  $T_3$  is bound to proteins (TBG and albumin). Abnormal levels (high or low) of thyroid hormone-binding proteins (primarily albumin and TBG) may cause abnormal  $T_3$  concentrations in euthyroid patients. This test is a measurement of total  $T_3$  (i.e., the free and the bound  $T_3$ ). Generally, when the  $T_3$  level is below normal, the patient is in a hypothyroid state.

In severe nonthyroidal diseases,  $T_3$  levels are decreased because conversion of  $T_4$  to  $T_3$  in the liver is diminished. This makes  $T_3$  levels less useful in indicating hypothyroid states. Instead, there is peripheral conversion of  $T_4$  to  $T_3$ . This form of  $T_3$  is called *reverse*  $T_3$  ( $rT_3$ ). Elevated  $rT_3$  levels associated with reduced  $T_3$  have been observed in starvation, anorexia nervosa, severe trauma and hemorrhagic shock, hepatic dysfunction, postoperative states, severe infection, and burn patients (i.e., “sick euthyroid” syndrome).

Furthermore, there is considerable overlap between hypothyroid states and normal thyroid function. Because of this,  $T_3$  levels are mostly used just to assist in the diagnosis of hyperthyroid states.  $T_3$  levels are frequently low in sick or hospitalized euthyroid patients.

An elevated  $T_3$  level indicates hyperthyroidism, especially when the  $T_4$  is also elevated. There is a rare form of hyperthyroidism called  $T_3$  *toxicosis*, in which the  $T_4$  is normal and the  $T_3$  is elevated.

## Interfering factors

- $T_3$  values are increased in pregnancy.
- ▼ Drugs that may cause *increased* levels include estrogen, methadone, and oral contraceptives.
- ▼ Drugs that may cause *decreased* levels include anabolic steroids, androgens, phenytoin, propranolol, reserpine, and salicylates (high dose).

## Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Determine whether the patient is taking any exogenous  $T_3$  medication because this will affect test results.
- Withhold drugs that may affect results (with physician's approval).

## Abnormal findings

### ▲ Increased levels

Graves disease  
 Plummer disease  
 Toxic thyroid adenoma  
 Acute thyroiditis  
 Factitious hyperthyroidism  
 Struma ovarii  
 Pregnancy  
 Hepatitis  
 Congenital  
 hyperproteinemia

### ▼ Decreased levels

Hypothyroidism  
 Congenital hypothyroidism  
 Thyroid surgical ablation  
 Myxedema  
 Pituitary insufficiency  
 Hypothalamic failure  
 Protein malnutrition and  
 other protein-depleted  
 states (e.g., nephrotic  
 syndrome)  
 Iodine insufficiency  
 Renal failure  
 Cushing syndrome  
 Cirrhosis  
 Liver diseases

notes

**troponins** (Cardiac-specific troponin T [cTnT], Cardiac-specific troponin I [cTnI], High-sensitivity troponin [hsTnT, hs-cTn])

**Type of test** Blood

### Normal findings

Cardiac troponin T: < 0.1 ng/mL

Cardiac troponin I: < 0.03 ng/mL

hsTnT :

< 14 ng/L for women

< 22 ng/L for men

### Test explanation and related physiology

Cardiac troponins are important biochemical markers for cardiac disease. This test is used to assist in the evaluation of patients with suspected acute coronary ischemic syndromes. In addition to improving the diagnosis of acute ischemic disorders, troponins are also valuable for early risk stratification in patients with unstable angina. They can be used to predict the likelihood of future cardiac events. There are two cardiac-specific troponins: cTnT and cTnI.

Because of their extraordinarily high specificity for myocardial cell injury, cardiac troponins are very helpful in the evaluation of patients with chest pain. Their use is similar to that of creatine kinase-MB (CK-MB). However, cardiac troponins have several advantages over CK-MB. Cardiac troponins are more specific for cardiac muscle injury. CK-MB can be elevated in severe skeletal muscle injury, brain or lung injury, or renal failure. Cardiac troponins will nearly always be normal in noncardiac muscle diseases. Cardiac troponins become elevated sooner and remain elevated longer than CK-MB. This expands the time window of opportunity for diagnosis and thrombolytic treatment of myocardial injury. Finally, cTnT and cTnI are more sensitive to muscle injury than CK-MB. That is most important in evaluating patients with chest pain.

Cardiac troponins become elevated as early as 2 to 3 hours after myocardial injury. Typically, 2 to 3 sets of troponins over the course of a day (every 3-6 hours) are required to indicate myocardial infarction. Levels of cTnI may remain elevated for 7 to 10 days after myocardial infarction, and cTnT levels may remain elevated for up to 14 days. High sensitivity (hs)TnT represents an important advance with added sensitivity for cardiac myocyte necrosis. hsTnT assays detect concentrations of troponins at much lower concentrations, thereby detecting infarction

earlier in the process (as early as 90 minutes). hsTnT speeds the triage of patients with suspected myocardial infarction.

Cardiac troponins are used in the following cardiac clinical situations:

- Evaluation of patient with unstable angina. If cardiac troponins are normal, no myocardial injury has occurred, and there will be no lasting cardiac dysfunction. If cardiac troponins are elevated, muscle injury has occurred. Revascularization may be indicated because this latter group is at great risk for a subsequent cardiac event (infarction or sudden death).
- Detection of reperfusion associated with coronary recanalization. A *washout*, or second peak of cardiac troponin levels, accurately indicates reperfusion by way of recanalization or coronary angioplasty.
- Estimation of myocardial infarction size. Late (4 weeks) cardiac troponin levels are inversely related to left ventricular ejection fraction.
- Detection of perioperative myocardial infarction. Cardiac troponins are not affected by skeletal muscle injury.
- Evaluation of the severity of pulmonary emboli. Elevated levels may indicate more severe disease and the need for thrombolytic therapy.
- Congestive heart failure. Persistently elevated troponins indicate continued ventricular strain.

Elevations of troponin T do not necessarily indicate the presence of an ischemic mechanism. Many other disease states are associated with elevations of troponin T via mechanisms different from those that cause injury in patients with acute coronary syndromes. These include cardiac trauma (e.g., contusion, ablation, or pacing), congestive heart failure, hypertension, hypotension (often with arrhythmias), pulmonary embolism, renal failure, and myocarditis.

### Interfering factors

- Severe skeletal muscle injury may cause false elevation of cTnT.
- Troponin T levels are falsely elevated in dialysis patients.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: yellow

**PT** Discuss with the patient the need and reason for frequent venipuncture in diagnosing myocardial infarction.

## 914 troponins

- After the initial blood sample, blood is collected 12 hours later followed by daily testing for 3 to 5 days and possibly weekly for 5 to 6 weeks.
- Record the exact time and date of venipuncture on each laboratory slip. This aids in the interpretation of the temporal pattern of enzyme elevations.
- If a qualitative immunoassay is to be done at the bedside, whole blood is obtained in a micropipette and placed in the sample well of the testing device. A red or purple color in the read zone indicates that 0.2 ng/mL or more of cardiac troponin is present in the patient's blood.

### Abnormal findings

#### ▲ Increased levels

Myocardial injury

Myocardial infarction

#### notes

## tuberculin skin test (TST, Tuberculin test, Mantoux test, Purified protein derivative [PPD] test)

**Type of test** Skin

### Normal findings

PPD: Negative; reaction < 5 mm

### Test explanation and related physiology

Tuberculin testing is performed for persons who are

- Suspected of having active TB
- At increased risk for progression to active TB
- At increased risk for latent TB infection (LTBI) (e.g., health-care workers, recent immigrants, or IV drug abusers)
- At low risk for LTBI, but are tested for other reasons (e.g., entrance to college)

For this test, a PPD of the tubercle bacillus is injected intradermally. If the patient is infected with TB (whether active or dormant), lymphocytes will recognize the PPD antigen and cause a local reaction; if the patient is not infected, no reaction will occur. If the test result is negative and the physician still strongly suspects TB, a *second-strength PPD* can be used. If this test result is negative, the patient does not have TB. (See p. 917 for TB culture.) A positive reaction usually occurs 6 weeks after infection. If positive, a localized skin reaction will recur with any subsequent testing throughout the person's life. This test is used to detect TB infection but is unable to indicate whether the infection is active or latent.

When a patient known to have active TB receives a PPD test, the local reaction may be so severe that it causes complete skin slough and requires surgical care. When these patients are eliminated from PPD testing, the test has no complications. The PPD test will not cause active TB because no live organisms exist in the test solution.

The PPD test also can be used as part of a series of skin tests to assess the immune system. If the immune system is nonfunctioning because of poor nutrition or chronic illness (e.g., neoplasia, infection), the PPD test will be negative despite the patient having had an active or dormant TB infection.

There are several alternative tests (see p. 919) for diagnosing *Mycobacterium tuberculosis* infection.

### Contraindications

- Patients with known active TB

## 916 tuberculin skin test

- Patients who have received *bacille Calmette-Guérin* (BCG) immunization against PPD, because these patients will demonstrate a positive reaction to the PPD vaccination even though they have never had TB infection

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Assure the patient that TB cannot result from testing.
  - Assess the patient for a previous history of TB. Report a positive history to the physician.
  - Evaluate the patient's history for previous PPD results and BCG immunization.

#### During

- Prepare the patient's forearm with alcohol and allow it to dry.
- Intradermally inject the PPD. A skin wheal will occur.
- Circle the area with indelible ink.
- Record the time at which the PPD was injected.

#### After

- Read the results in 48 to 72 hours.
- Examine the test site for induration (hardening). Measure the area of induration (not redness) in millimeters.
- If the test is positive, ensure that the physician is notified and the patient is treated appropriately.
- If the test is positive, check the patient's arm 4 to 5 days after the test to be certain that a severe skin reaction has not occurred.

### Abnormal findings

#### Positive results

- TB, infection
- Nontuberculous mycobacteria infections

#### Negative results

- Possible immunoincompetence

#### notes



**tuberculosis culture** (TB culture, Acid-fast bacillus [AFB] smear)

**Type of test** Microbiology culture

**Normal findings**

Negative for TB

**Test explanation and related physiology**

Diagnosis of TB can be made only by identification and culture of *Mycobacterium tuberculosis* in a specimen. Conventional culture techniques for growth, identification, and susceptibility of acid-fast mycobacteria take 4 to 6 weeks. Because a patient suspected of having TB cannot be isolated for that duration, the disease may spread to many other people while he or she waits for the diagnosis. With the resurgence and increasing incidence of TB in the U.S. population (especially among immunocompromised patients with AIDS), newer, more rapid culture techniques have been identified and are now being performed.

The *BACTEC* method is a culture technique in which the growth medium for culturing mycobacteria is supplanted with a substrate labeled with radioactive carbon ( $^{14}\text{C}$ ). This substrate is used by mycobacteria; during metabolism, radioactive carbon dioxide ( $^{14}\text{CO}_2$ ) is produced from the substrate. The  $^{14}\text{CO}_2$  is detected and quantitated. This permits quick identification of mycobacterial growth.

*PCR culture* methods have also been developed. With the addition of a DNA polymerase, genetic chromosomal parts can be multiplied. This allows amplification of genomes, which then can be detected by genetic DNA probes. With these newer techniques, *M. tuberculosis* can be identified in as few as 36 to 48 hours. With this reduction in diagnostic time, treatment can be started earlier.

When TB is suspected, a sputum smear for *AFB* can be obtained. After taking up a dye, such as fuchsin, *M. tuberculosis* is not decolorized by acid alcohol (i.e., it is acid fast). It is seen under the microscope as a red, rod-shaped organism. If this bacillus is seen, the patient may have active TB. Other specimens, such as cerebrospinal fluid, tissue, and synovial fluid, may be used. *AFB* is also used to monitor treatment for TB. If after adequate therapy (2 months), the sputum still contains *AFB* (even though the culture may be negative because of anti-TB drugs), treatment failure should be considered.

## Interfering factors

- ▼ Antituberculosis drugs

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Tell the patient that no fasting is required.

### During

- For *sputum* collection, it is best to induce sputum production with an ultrasonic or nebulizing device. Collect three to five early morning sputum specimens. All specimens must contain mycobacteria to make the diagnosis of TB.
- For *urine* collection, obtain three to five single clean-voided specimens early in the morning.
- *Swabs, intestinal washings, and biopsy* specimens should be transported to the laboratory immediately for preparation.
- When the specimen is received by the laboratory, a decontamination process is applied to it to kill all nonmycobacteria. The specimen is then cultured in the appropriate medium.
- With the rapid growth techniques, the specimen is evaluated every 24 hours.
- When culture growth is considered adequate, the organisms are stained for AFB and identified.
- With genetic probes, the *Mycobacterium* species is identified.
- At this point, if *M. tuberculosis* is present, the report will read “culture more positive for mycobacteria.” If the species has been identified, this also will be reported.
- Drug-susceptibility testing then will be carried out and subsequently reported.

### After

- PT** Instruct the patient as to appropriate isolation of sputum and other body fluids to avoid potential spread of suspected TB.

## Abnormal findings

Tuberculosis

Atypical mycobacterial nontuberculosis disease

notes

**tuberculosis testing** (TB testing, Interferon gamma release assay [IGRA, T-Spot.TB], QuantiFERON-TB Gold [QFT, QFT-G, TB Gold test, TB blood test], Nucleic acid amplification for TB [NAAT], TB antibody)

**Type of test** Blood, sputum, fluid analysis

### Normal findings

IGRA result	Interpretation
Positive	<i>Mycobacterium tuberculosis</i> infection likely
Negative	<i>M. tuberculosis</i> infection unlikely but cannot be excluded. If TB disease is highly suspected, a negative result does not rule out infection. False-negative results may be seen in immunocompromised patients.
Indeterminate	Test not interpretable. Collection of a new specimen for testing is recommended, if clinically indicated.

### Test explanation and related physiology

Blood tests such as interferon gamma release assays (IGRA, e.g., QuantiFERON-TB or T-Spot.TB), nucleic acid amplification test (NAAT), and serologic TB testing are used to diagnose active or latent TB infection. The gold standard for making the diagnosis of active TB is the TB culture (see p. 917). However, it takes 2 to 6 weeks to obtain results. Identifying acid-fast bacilli in a smear (AFB smear) (see p. 917) of the body fluid (usually sputum) is a rapid method of identifying TB in 24 hours. Unfortunately, AFB is not very sensitive or specific. The IGRA is a whole-blood test used in diagnosing *Mycobacterium tuberculosis* infection. The NAAT is a rapid and accurate test of sputum and is used as corroborative information in the diagnosis of TB. Serology testing on blood is also a rapid test used to identify active TB disease infection.

The determination as to whether the patient has active or latent TB still requires additional testing (including chest radiography, sputum smear, and culture). IGRA is a diagnostic aid that measures a component of cell-mediated immune reactivity to *M. tuberculosis* similar to the tuberculin skin testing (TST) (see p. 915). IGRA can be performed on patients with prior *bacille Calmette-Guérin* (BCG) vaccination without causing a hypersensitivity response. Similar to TST, false negatives can occur in anergic patients.

## 920 tuberculosis testing

IGRA and NAAT are used in the same patient population as TST. These include contact investigations, evaluation of recent immigrants (from some countries), HIV, dialysis, malnourished patients, and sequential-testing surveillance programs for infection control, such as those for health-care workers.

NAAT is designed to identify TB complex DNA in a body fluid (bronchoalveolar lavage, bronchial washing, sputum, stool, pleural or abdominal fluid, tissue, or urine sample). This test provides a rapid result in 24 hours. Not only does this test provide early diagnosis and allow treatment, it can also indicate resistance to rifampicin, a drug commonly used to treat TB. Similar to the above-described rapid tests, NAAT cannot indicate active infection from a previously treated TB infection.

TB antibody serology is designed to identify IgG antibodies to TB mycobacteria in patients with active TB infections. This blood test can be used in previously vaccinated BCG patients. It is particularly useful in evaluating the effectiveness of anti-TB therapy and documenting a response to therapy. Similar to IGRA, serology results may not be positive in immunocompromised patients, making use in HIV-infected patients less helpful.

### Interfering factors

- Heterophile (e.g., human antimouse) antibodies in serum or plasma of certain individuals are known to cause interference with immunoassays.
- A false-negative IGRA result can be because of the stage of infection (i.e., specimen obtained before the development of cellular immune response), comorbid conditions that affect immune function, or other individual immunologic factors.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient or the family.

#### During

- Collect 1 mL of whole blood in each of three laboratory-specified collection tubes. The accuracy of the IGRA depends on the proper collection and incubation of the blood specimen. Blood should fill the tube as close to the 1-mL mark as possible. Under- or overfilling the tubes outside the 0.8- to 1.2-mL range may lead to erroneous results.
- Immediately after collection, each specimen tube must be shaken vigorously by shaking the tube up and down 10 times to ensure that the entire inner surface of the tube has been coated with blood. This distributes the stimulating antigens,

allowing optimal processing and presentation of the antigens to T-cells, which causes release of interferon- $\gamma$ .

- For NAAT testing, 1 to 3 mL of sputum or body fluid is required. This should be refrigerated in a screw cap sterile container.

### After

- Apply pressure or a pressure dressing to the venipuncture site.
- PT** If the patient's results are positive, educate the patient on the necessary follow-up studies, such as chest radiography and sputum cultures.

### Abnormal finding

TB infection

notes

**upper gastrointestinal x-ray study** (Upper GI series, UGI)

**Type of test** X-ray with contrast dye

**Normal findings**

Normal size, contour, patency, filling, positioning, and transit of barium through the lower esophagus, stomach, and upper duodenum

**Test explanation and related physiology**

The upper GI study consists of a series of x-ray images of the lower esophagus, stomach, and duodenum, usually with barium sulfate as the contrast medium. When there is concern about leakage of x-ray contrast through a perforation of the GI tract, however, Gastrografin (a water-soluble contrast) is used. This test can be performed in conjunction with a barium swallow or a small bowel series (pp. 131 and 830), which can precede or succeed the upper GI study, respectively.

The purpose of this examination is to detect ulcerations, tumors, inflammations, and anatomic malpositions (e.g., hiatal hernia) in these organs. Obstruction of the upper GI tract is also easily detected.

In this test, the patient is asked to drink barium. As the contrast descends, the lower esophagus is examined for position, patency, and filling defects (e.g., tumors, scarring, varices). As the contrast enters the stomach, the gastric wall is examined for benign or malignant ulcerations, filling defects (most often in cancer), and anatomic abnormalities (e.g., hiatal hernia). The patient is placed in a flat or head-down position, and the gastroesophageal area is examined for evidence of gastroesophageal reflux of barium.

As the contrast leaves the stomach, patency of the pyloric channel and the duodenum is evaluated. Benign peptic ulceration is the most common pathologic condition affecting these areas. Extrinsic compression caused by tumors, cysts, or enlarged pathologic organs (e.g., liver) near the stomach also can be identified by anatomic distortion of the outline of the upper GI tract. The small intestine can then be studied (see discussion of *small bowel follow-through*, p. 830).

**Contraindications**

- Patients with complete bowel obstructions
- Patients suspected of upper GI perforation. Water-soluble Gastrografin should be used instead of barium.

- Patients with unstable vital signs. These patients should be supervised during testing.
- Patients who are uncooperative because of the necessity of frequent position changes

### Potential complications

- Aspiration of barium
- Constipation or partial bowel obstruction caused by inspissated barium in the small bowel or colon

### Interfering factors

- Previously administered barium may block visualization.
- Incapacitated patients cannot assume the multiple positions required for the study.
- Food and fluid in the stomach give the false impression of filling defects within the stomach, precluding adequate evaluation of the gastric mucosa.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. Allow the patient to verbalize concerns. See p. xviii for radiation exposure and risks.
- PT** Instruct the patient to abstain from eating for at least 8 hours before the test. Usually keep the patient NPO (nothing by mouth) after midnight on the day of the test.
- PT** Assure the patient that the test will not cause any discomfort.

#### During

- Note the following procedural steps:
  1. The patient is asked to drink approximately 16 oz of barium sulfate. This is a chalky substance usually suspended in milk shake form and drunk through a straw. Usually the drink is flavored to increase palatability.
  2. After drinking the barium, the patient is moved through several position changes (e.g., prone, supine, lateral) to promote filling of the entire upper GI tract.
  3. Images are taken at the discretion of the radiologist observing the flow of barium fluoroscopically.
  4. The flow of barium is followed through the lower esophagus, stomach, and duodenum.
  5. In an *air-contrast upper GI study*, the patient is asked to rapidly swallow carbonated powder. This creates carbon dioxide in the stomach, providing air contrast to the barium in the stomach and increased visualization of the gastric mucosa.

## 924 upper gastrointestinal x-ray study

- Note that a radiologist performs this procedure in approximately 30 minutes.
- PT** Tell the patient that he or she may be uncomfortable lying on the hard x-ray table and may experience the sensation of bloating or nausea during the test.

### After

- PT** Inform the patient that if Gastrografin was used, he or she may have significant diarrhea. Gastrografin is an osmotic cathartic.
- PT** Instruct the patient to use a cathartic (e.g., milk of magnesia) if barium sulfate was used as the contrast medium. Water absorption may cause the barium to harden and create fecal impaction if catharsis is not carried out.
- PT** Instruct the patient to watch his or her stools to ensure that all of the barium has been removed. The stools should return to a normal color after completely expelling the barium, which may take as long as a day and a half.

### Abnormal findings

Esophageal cancer

Esophageal varices

Hiatal hernia

Esophageal diverticula

Gastric cancer

Gastric inflammatory disease (e.g., Ménétrier disease)

Benign gastric tumor (e.g., leiomyoma)

Extrinsic compression by pancreatic pseudocysts, cysts, pancreatic tumors, or hepatomegaly

Perforation of the esophagus, stomach, or duodenum

Congenital abnormalities (e.g., duodenal web, pancreatic rest, malrotation syndrome)

Gastric ulcer (benign and malignant)

Gastritis

Duodenal ulcer

Duodenal cancer

Duodenal diverticulum

### notes



## urea breath test (UBT, *Helicobacter pylori* breath test)

**Type of test** Other

### Normal findings

< 50 dpm (if  $^{14}\text{C}$  is used)

< 3% (if  $^{13}\text{C}$  is used)

### Test explanation and related physiology

This test is used to detect *Helicobacter pylori* infections. It is indicated in patients who have recurrent or chronic gastric or duodenal ulceration or inflammation. When the *H. pylori* infection is successfully treated, the ulcer or inflammation will usually heal.

There are several serologic and microscopic methods of detecting *H. pylori* (see *Helicobacter pylori* antibodies test, p. 482). The UBT is the noninvasive test of choice for diagnosis of *H. pylori* infection. It is based on the capability of *H. pylori* to metabolize urea to  $\text{CO}_2$  because of the organism's capability to produce a large amount of urease. In the breath test, carbon ( $^{13}\text{C}$ )-labeled urea is administered orally. The urea is then absorbed through the gastric mucosa. If *H. pylori* is present, the urea will be converted to  $^{13}\text{CO}_2$ . The  $^{13}\text{CO}_2$  is then taken up by the capillaries in the stomach wall and delivered to the lungs, where it is exhaled. The labeled carbon can be measured by gas chromatography or a mass spectrometer.

This test has been simplified to the point that two breath samples collected before and 30 minutes after the ingestion of urea in a liquid form suffice to provide reliable diagnostic information. Labeling urea with  $^{13}\text{C}$  is becoming increasingly popular because it is a nonradioactive isotope of  $^{14}\text{C}$  and is innocuous. It can be safely used in children and women of childbearing age.

### Interfering factors

- Dietary constituents with a natural abundance of  $^{13}\text{C}$ , such as maize, cane, and corn flour, can cause increased levels.
- ✚ Bismuth (Pepto-Bismol) or sucralfate (Carafate) will suppress mucosal uptake of the urea and interfere with test results.
- ✚ The concomitant use of a proton pump inhibitor, such as Prilosec, Nexium, Prevacid, or Protonix, will also inhibit urea absorption.

## 926 urea breath test

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the patient to abstain from oral intake for 6 hours before testing.
  - If radioactive carbon (rare) is being used, be sure that female patients are not pregnant.
- PT** When providing the isotopic urea to the patient, instruct the patient as to proper administration (per laboratory routine).

#### During

- Several minutes after the patient has swallowed the carbon dose, provide the patient with 2 oz of water.
- Breath samples are collected in any one of a number of gas collection devices depending on how and when the sample will be analyzed.

#### After

- PT** Instruct the patient to resume medications and normal diet.
- PT** If radioactive carbon was used, instruct the patient to drink plenty of fluids to facilitate excretion of the radioisotope.

### Abnormal findings

*H. pylori* infection

notes

## uric acid, blood, and urine

**Type of test** Blood, urine

### Normal findings

#### Blood

##### Adult

Male: 4.0-8.5 mg/dL or 0.24-0.51 mmol/L

Female: 2.7-7.3 mg/dL or 0.16-0.43 mmol/L

Elderly: values may be slightly increased

Child: 2.5-5.5 mg/dL or 0.12-0.32 mmol/L

Newborn: 2.0-6.2 mg/dL

Physiologic saturation threshold: > 6 mg/dL or > 0.357 mmol/L

Therapeutic target for gout: < 6 mg/dL or < 0.357 mmol/L

#### Urine

250-750 mg/24 hr or 1.48-4.43 mmol/day (SI units)

### Possible critical values

Blood: > 12 mg/dL

### Test explanation and related physiology

Uric acid is a nitrogenous compound that is a product of purine (a DNA building block) catabolism. Uric acid is excreted to a large degree by the kidney and to a smaller degree by the intestinal tract. When uric acid levels are elevated (hyperuricemia), the patient may have gout. Gout is a common metabolic disorder characterized by chronic hyperuricemia. At high levels, uric acid concentrations exceed the physiologic saturation threshold, and monosodium urate crystals may be deposited in the joints and soft tissues.

Causes of hyperuricemia can be overproduction or decreased excretion of uric acid (e.g., kidney failure). Overproduction of uric acid may occur in patients with a catabolic enzyme deficiency that stimulates purine metabolism or in patients with cancer in whom purine and DNA turnover is great. Other causes of hyperuricemia may include alcoholism, leukemia, metastatic cancer, multiple myeloma, hyperlipoproteinemia, diabetes mellitus, renal failure, stress, lead poisoning, and dehydration caused by diuretic therapy. Ketoacids (as occur in diabetic or alcoholic ketoacidosis) may compete with uric acid for tubular excretion and may cause decreased uric acid excretion. Many causes of hyperuricemia are undefined and therefore labeled as *idiopathic*.

## 928 uric acid, blood, and urine

Elevated uric acid in the urine is called uricosuria. Uric acid can become supersaturated in the urine and crystallize to form kidney stones that can block the renal system. Urinary excretion of uric acid depends on uric acid levels in the blood, along with glomerular filtration and tubular secretion of uric acid into the urine. Uric acid is less well saturated in alkaline urine. As the urine pH rises, more uric acid can exist without crystallization and stone formation. Therefore when a person is known to have high uric acid in the urine, the urine can be alkalinized by ingestion of a strong base to prevent stone formation.

### Interfering factors

- Stress may cause increased uric acid levels.
- Recent use of x-ray contrast agents may cause decreased serum levels.
- Recent use of x-ray contrast agents may increase uric acid levels in the urine.
- ✔ Drugs that may cause *increased* serum levels include alcohol, ascorbic acid, aspirin (low dose), caffeine, cisplatin, diazoxide, diuretics, epinephrine, ethambutol, levodopa, methyl dopa, nicotinic acid, phenothiazines, and theophylline.
- ✔ Drugs that may cause *decreased* serum levels include allopurinol, aspirin (high dose), azathioprine, clofibrate, corticosteroids, estrogens, glucose infusions, guaifenesin, mannitol, probenecid, and warfarin.
- ✔ Drugs that may cause *increased* urine levels include ascorbic acid, calcitonin, citrate, dicumarol, estrogens, glyceryl, iodinated dyes, phenolsulfonphthalein, probenecid, salicylates, steroids, and outdated tetracycline.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- Follow the institution's requirements regarding fasting.

#### During

##### *Blood*

- Collect a venous blood sample in a red-top tube.

##### *Urine*

- PT** Instruct the patient to begin the 24-hour urine collection after voiding. Discard the initial specimen and start the 24-hour timing at that point.
- See inside front cover for Routine Urine Testing.

**After***Blood*

- Apply pressure to the venipuncture site.

*Urine*

- Transport the urine specimen promptly to the laboratory.

**Abnormal findings****▲ Increased blood levels (hyperuricemia)**

Gout  
 Increased ingestion of purines  
 Genetic inborn error in purine metabolism  
 Metastatic cancer  
 Multiple myeloma  
 Leukemia  
 Cancer chemotherapy  
 Hemolysis  
 Rhabdomyolysis (e.g., heavy exercise, burns, crush injury, epileptic seizure, or myocardial infarction)  
 Chronic renal disease  
 Acidosis (ketotic or lactic)  
 Hypothyroidism  
 Toxemia of pregnancy  
 Hyperlipoproteinemia  
 Alcoholism  
 Shock or chronic blood volume depletion states  
 Idiopathic

**▲ Increased urine levels**

Gout  
 Metastatic cancer  
 Multiple myeloma  
 Leukemia  
 Cancer chemotherapy  
 High purine diet  
 Lead toxicity

**▼ Decreased blood levels**

Wilson disease  
 Fanconi syndrome  
 Lead poisoning  
 Yellow atrophy of the liver

**▼ Decreased urine levels**

Kidney disease  
 Eclampsia  
 Chronic alcohol ingestion  
 Acidosis (ketotic or lactic)

**urinalysis (UA)****Type of test** Urine**Normal findings**

Appearance: clear

Color: amber yellow

Odor: aromatic

pH: 4.6 to 8 (average 6)

Protein

0-8 mg/dL

50-80 mg/24 hr (at rest)

&lt; 250 mg/24 hr (exercise)

Specific gravity

Adult: 1.005-1.030 (usually 1.010-1.025)

Elderly: values decrease with age

Newborn: 1.001-1.020

Leukocyte esterase: negative

Nitrites: negative

Ketones: negative

Crystals: negative

Casts: none present

Glucose

Fresh specimen: negative

24-hour specimen: 50-300 mg/day or 0.3-1.7 mmol/day (SI units)

White blood cells (WBCs): 0-4 per low-power field

WBC casts: negative

Red blood cells (RBCs):  $\leq 2$ 

RBC casts: none

**Test explanation and related physiology**

A total urinalysis involves multiple routine tests on a urine specimen. This specimen is not necessarily a clean-catch specimen. However, if urinary tract infection (UTI) is suspected, often a midstream clean-catch specimen is obtained. This urine is then split into two parts. One is sent for urinalysis, and the other is held in the laboratory refrigerator and cultured (see p. 947) if the urinalysis indicates infection. Abnormalities detected by urinalysis may reflect either urinary tract diseases (e.g., infection, glomerulonephritis, loss of concentrating capacity) or extrarenal disease processes (e.g., glucosuria in diabetes, proteinuria in monoclonal gammopathies, bilirubinuria in liver disease).

Urinalysis routinely includes remarks about the color, appearance, and odor of the urine. The pH is determined. The urine is tested for the presence of proteins, glucose, ketones, blood, and leukocyte esterase. The urine is examined microscopically for RBCs, WBCs, casts, crystals, and bacteria.

Examination of the urine sediment provides a significant amount of information about the urinary system. Reference ranges have been provided to recognize marked abnormalities.

### Appearance and color

Urine appearance and color are noted as part of routine urinalysis. The appearance of a normal urine specimen should be clear. Cloudy urine may be caused by the presence of pus, RBCs, or bacteria; however, normal urine also may be cloudy because of ingestion of certain foods (e.g., large amounts of fat, urates, or phosphates). The color of urine ranges from pale yellow to amber because of the pigment *urochrome*. The color indicates the concentration of the urine and varies with specific gravity. Dilute urine is straw colored, and concentrated urine is deep amber.

Abnormally colored urine may result from a pathologic condition or the ingestion of certain foods or medicines. For example, bleeding from the kidney produces dark red urine, whereas bleeding from the lower urinary tract produces bright red urine. Dark yellow urine may indicate the presence of urobilinogen or bilirubin. *Pseudomonas* organisms may produce green urine. Beets may cause red urine, and rhubarb can color the urine brown. Many frequently used drugs also may affect urine color (Table 37).

### Odor

Determination of urine odor is a part of routine urinalysis. The aromatic odor of fresh, normal urine is caused by the presence of volatile acids. Urine of patients with diabetic ketoacidosis has the strong, sweet smell of acetone. In patients with UTIs, the urine may have a very foul odor. Patients with a fecal odor to their urine may have an enterovesical fistula.

### pH

Urine pH is affected by diet, medications, systemic acid-base disturbances, and renal tubular function. An alkaline pH is obtained in a patient with alkalemia. Also, bacteria, UTI, or a diet high in citrus fruits or vegetables may cause an increased urine pH. An alkaline urine pH is common after eating. Certain medications (e.g., streptomycin, neomycin, kanamycin) are effective in treating UTIs when the urine is alkaline. Slightly

**TABLE 37** Frequently used drugs that may affect urine color

Generic and brand names	Drug classification	Urine color
Cascara sagrada	Stimulant laxative	Red in alkaline urine; yellow-brown in acid urine
Chloroquine (Aralen)	Antimalarial	Rusty yellow or brown
Chlorzoxazone (Paraflex)	Skeletal muscle relaxant	Orange or purple-red
Docusate calcium (Doxidan, Surfak)	Laxative	Pink to red to red-brown
Doxorubicin (Adriamycin)	Antineoplastic	Red-orange
Iron preparations (Ferrotran, Imferon)	Hematinic	Dark brown or black on standing
Levodopa	Antiparkinsonian	Dark brown on standing
Metronidazole (Flagyl)	Antiinfective	Darkening, red-brown
Nitrofurantoin (Macrochantin, Nitrofan)	Antibacterial	Brown-yellow
Phenazopyridine (Pyridium)	Urinary tract analgesic	Orange to red
Phenolphthalein (Ex-Lax)	Contact laxative	Red or purple-pink in alkaline urine
Phenothiazines (e.g., prochlorperazine [Compazine])	Antipsychotic, neuroleptic, antiemetic	Red-brown
Phenytoin (Dilantin)	Anticonvulsant	Pink, red, red-brown
Riboflavin (vitamin B)	Vitamin	Intense yellow
Rifampin	Antibiotic	Red-orange
Sulfasalazine (Azulfidine)	Antibacterial	Orange-yellow in alkaline urine
Triamterene (Dyrenium)	Diuretic	Pale blue fluorescence



acidic urine is normal. Normal average pH is 7.0, slightly acidic compared with average blood pH of 7.4. However, acidic urine is observed in patients with acidemia, which can result from metabolic or respiratory acidosis, starvation, dehydration, or a diet high in meat products or cranberries. In patients with renal tubular acidosis, however, the blood is acidic, and the urine is alkaline.

The urine pH is useful in identifying crystals in the urine and determining the predisposition to form a given type of stone. Acidic urine is associated with xanthine, cystine, uric acid, and calcium oxalate stones. To treat or prevent these urinary calculi, urine should be kept alkaline. Alkaline urine is associated with calcium carbonate, calcium phosphate, and magnesium phosphate stones; for these stones urine should be kept acidic. (See *stone analysis*, p. 945.)

### Protein

Protein is a sensitive indicator of glomerular and tubular renal function. Normally, less than 30 mg of protein per day is in urine. This is not detectable in the routine protein analysis. Microalbumin can be detected, however (see p. 626). If the glomerular filtration membrane is injured, as in glomerulonephritis, the spaces become much larger, and protein seeps out into the filtrate and then into the urine. Renal tubules are a site of protein reabsorption. If tubular disease exists, protein is in the urine. If proteinuria persists at a significant rate, the patient can become hypoproteinemic because of the severe protein loss through the kidneys. This decreases the normal capillary oncotic pressure that holds fluid in the vasculature and causes severe interstitial edema. The combination of proteinuria and edema is known as the nephrotic syndrome.

Proteinuria (usually albumin) is probably the most important indicator of renal disease. The urine of all pregnant women is routinely checked for proteinuria, which can be an indicator of preeclampsia. In addition to screening for nephrotic syndrome, urinary protein also screens for complications of diabetes mellitus, glomerulonephritis, amyloidosis, and multiple myeloma (see test for Bence Jones protein, p. 134).

If significant protein is noted at urinalysis, a 24-hour urine specimen should be collected to measure the quantity of protein. This estimate of 24-hour protein excretion is usually performed with a urine creatinine because hydration status and other factors may influence urine concentration. The normal *protein/creatinine ratio* is less than 0.15.

**Glucose**

This can be an effective screening test for the presence of glucose in the urine that may identify diabetes mellitus or other causes of glucose intolerance (see glucose, p. 467). Although urine glucose tests previously were used to monitor the effectiveness of diabetes therapy, today glucose monitoring is largely done by fingerstick determinations of blood glucose levels.

Glucose is filtered from the blood by the glomeruli of the kidney. Normally, all of the glucose is resorbed in the proximal renal tubules. When the blood glucose level exceeds the capability of the renal threshold to resorb the glucose (normally ~180 mg/dL), it begins to spill over into the urine (glycosuria). As the blood glucose level increases further, greater amounts of glucose are spilled into the urine.

Glucosuria may occur immediately after eating a high-carbohydrate meal in patients with a low tubular maximum for glucose. Similarly, glucosuria can occur with a normal serum glucose level when kidney disease affects the renal tubule. The renal threshold for glucose becomes abnormally low, and glucosuria occurs. Glucosuria is not abnormal, however, in patients receiving dextrose-containing intravenous fluids. Patients with acute severe physical stress or injury can have a transient glucosuria caused by normal compensatory endocrine-mediated responses.

**Osmolality**

This is measured during routine urinalysis (see p. 656).

**Specific gravity**

Specific gravity is a measure of the concentration of particles, including wastes and electrolytes, in the urine. A high specific gravity indicates concentrated urine; a low specific gravity indicates dilute urine. Specific gravity refers to the weight of the urine compared with that of distilled water (which has a specific gravity of 1.000). Particles in the urine give it weight or specific gravity.

Specific gravity is used to evaluate the concentrating and excretory power of the kidney. Renal disease tends to diminish the concentrating capability of the kidney. As a result, chronic renal diseases are associated with a low specific gravity. Nephrogenic *diabetes insipidus* is associated with very little variation in specific gravity of the urine because the kidney cannot respond to such variables as hydration and solute load. Specific gravity is also a measurement of the hydration status of the patient. An overhydrated patient will have a more dilute urine with a lower specific gravity. The specific gravity of the urine in a dehydrated patient

can be expected to be abnormally high. Specific gravity correlates roughly with osmolality (see p. 656).

### **Leukocyte (WBC) esterase**

Leukocyte (WBC) esterase is a screening test used to detect leukocytes in the urine. When positive, this test indicates a UTI. This examination uses chemical testing with a leukocyte esterase dipstick; a shade of purple is considered a positive result. Some laboratories have established screening protocols in which a microscopic examination (see p. 947) is performed only if a leukocyte esterase test is positive.

### **Nitrites**

Like the leukocyte esterase test, the nitrite test is a screening test for the identification of UTIs. This test is based on the principle that many, but not all, bacteria produce an enzyme called *reductase*, which can reduce urinary nitrates to nitrites. Chemical testing is done with a dipstick containing a reagent that reacts with nitrites to produce a pink color, thus indirectly suggesting the presence of bacteria. A positive test result indicates the need for a urine culture. Nitrite screening enhances the sensitivity of the leukocyte esterase test to detect UTIs.

### **Ketones**

Normally no ketones are present in the urine; however, a patient with poorly controlled diabetes who is hyperglycemic may have massive fatty acid catabolism. Ketones (beta-hydroxybutyric acid, acetoacetic acid, and acetone) are the end products of this fatty acid breakdown. As with glucose, ketones (predominantly acetoacetic acid) spill over into the urine when the blood levels of patients with diabetes are elevated. The excess production of ketones in the urine is usually associated with poorly controlled diabetes. This test for ketonuria is also important in evaluating ketoacidosis associated with alcoholism, fasting, starvation, high-protein diets, and isopropanol ingestion. Ketonuria may occur with acute febrile illnesses, especially in infants and children.

### **Bilirubin and urobilinogen**

Bilirubin is a major constituent of bile. If bilirubin excretion is inhibited, conjugated (direct) hyperbilirubinemia will result (see p. 137). Unlike the unconjugated form, conjugated bilirubin is water soluble and can be excreted into the urine. Therefore bilirubin in urine suggests disease affecting bilirubin metabolism after conjugation or with defects in excretion (e.g., gallstones). Unconjugated bilirubin caused by prehepatic jaundice will not be excreted in the urine because it is not water soluble.

Bilirubin is excreted by way of the bile ducts into the bowel. There, some of the bilirubin is transformed into *urobilinogen* by the action of bacteria in the bowel. Most of the urobilinogen is excreted from the liver back into the bowel, but some is excreted by the kidneys.

### **Crystals**

Crystals found in urinary sediment on microscopic examination indicate that renal stone formulation is imminent, if not already present. Urea crystals occur in patients with high serum uric acid levels (gout). Phosphate and calcium oxalate crystals occur in the urine of patients with parathyroid abnormalities or malabsorption states. The type of crystal found varies with the disease and the pH of the urine.

### **Casts**

Casts are rectangular clumps of materials or cells that are formed in the renal distal and collecting tubules, where the material is maximally concentrated. These clumps of material and cells take on the shape of the tubule, thus the term *cast*. Casts are usually associated with some degree of proteinuria and stasis in the renal tubules. There are two kinds of casts: hyaline and cellular.

*Hyaline casts* are conglomerations of protein and indicate proteinuria. A few hyaline casts are normally found, especially after strenuous exercise or dehydration.

*Cellular casts*, which are conglomerations of degenerated cells, are described in the following paragraphs.

*Granular casts* result from the disintegration of cellular material into granular particles within a WBC or epithelial cell cast. Granular casts are found after exercise and in patients with various renal diseases.

In some diseases, the epithelial cells desquamate into the renal tubule. As the cells degenerate, fatty deposits in the cells coalesce and become incorporated with protein into *fatty casts*. These casts are all associated with glomerular disease or the nephrotic syndrome or nephrosis. Free oval fat bodies may also be associated with fatty emboli that occur in patients with bone fractures.

*Waxy casts* may be cell casts, hyaline casts, or renal failure casts. Waxy casts probably represent further degeneration of granular casts. They occur when urine flow through the renal tubule is diminished, giving time for granular casts to degenerate. Waxy casts are found especially in patients with chronic renal diseases and are associated with chronic renal failure. They also occur in patients with diabetic nephropathy, malignant hypertension, and glomerulonephritis.

Epithelial cells can enter the urine anywhere along the process of urinary excretion. The presence of occasional epithelial cells is not remarkable. Large numbers, however, are abnormal and can conglomerate into *tubular (epithelial) casts*. These are most suggestive of renal tubular disease or toxicity.

Normally, few WBCs are found in urine sediment on microscopic examination. The presence of five or more WBCs in the urine indicates a UTI involving the bladder, kidney, or both. A clean-catch urine culture should be done for further evaluation. *WBC casts* are most commonly found in infections of the kidney, such as acute pyelonephritis.

Any disruption in the blood–urine barrier, whether at the glomerular, tubular, or bladder level, will cause RBCs to enter the urine. The bleeding can be microscopic or gross hematuria. Patients with more than three RBCs per high-power field in two of three properly collected urine specimens should be considered to have microhematuria and hence be evaluated for possible pathologic causes. *RBC casts* suggest glomerulonephritis. RBC casts are also seen in patients with acute tubular necrosis, pyelonephritis, renal trauma, or renal tumor.

## Interfering factors

### Appearance and color

- Sperm in the urethra can cause the urine to appear cloudy.
- Urine that has been refrigerated for longer than 1 hour can become cloudy.
- Certain foods affect urine color: carrots may cause dark yellow urine; beets may cause red urine; rhubarb may cause reddish or brownish urine.
- Urine darkens with prolonged standing because of oxidation of bilirubin metabolites.
- ✚ Many drugs, given the right environment, can alter the color of urine (see [Table 37](#)).

### Odor

- Some foods (e.g., asparagus) produce a characteristic odor.
- When urine stands for a long time and begins to decompose, it has an ammonia-like smell.

### pH

- Urine pH becomes alkaline on standing, because of the action of urea-splitting bacteria, which produce ammonia.
- The urine pH of an uncovered specimen will become alkaline because carbon dioxide vaporizes from the urine.

## 938 urinalysis

- Dietary factors affect urine pH. Whereas ingestion of large quantities of citrus fruits, dairy products, and vegetables produces alkaline urine, a diet high in meat and certain foods (e.g., cranberries) produces acidic urine.
- ✚ Drugs that *increase* urine pH include acetazolamide, bicarbonate antacids, and carbonic anhydrase inhibitors.
- ✚ Drugs that *decrease* urine pH include ammonium chloride, chlorothiazide, and mandelic acid.

### Protein

- Transient proteinuria may be associated with severe emotional stress, excessive exercise, and cold baths.
- Radiopaque contrast media administered within 3 days may cause false-positive results for proteinuria.
- Urine contaminated with prostatic or vaginal secretions commonly causes proteinuria.
- Diets high in protein can cause proteinuria.
- Highly concentrated urine may have a higher concentration of protein than more dilute urine.
- Hemoglobin may cause a positive result with the dipstick method.
- Bence Jones protein may not appear with the dipstick method.
- ✚ Drugs that may cause *increased* protein levels include acetazolamide, aminoglycosides, amphotericin B, cephalosporins, colistin, griseofulvin, lithium, methicillin, nafcillin, nephrotoxic drugs, oxacillin, penicillamine, penicillin G, phenazopyridine, polymyxin B, salicylates, sulfonamides, tolbutamide, and vancomycin.

### Specific gravity

- Recent use of radiographic dyes increases specific gravity.
- Cold temperatures cause falsely high specific gravity.
- ✚ Drugs that may cause *increased* specific gravity include dextran, mannitol, and sucrose.

### Leukocyte esterase

- False-positive results may occur in specimens contaminated by vaginal secretions (e.g., heavy menstrual discharge, *Trichomonas* infection, parasites) that contain WBCs.
- False-negative results may occur in specimens containing high levels of protein or ascorbic acid.

### Ketones

- Special diets (carbohydrate-free, high-protein, high-fat) may cause ketonuria.

- ✚ Drugs that may cause false-positive results include bromosulfophthalein, isoniazid, isopropanol, levodopa, paraldehyde, phenazopyridine, and phenolsulfonphthalein.

### Bilirubin and urobilinogen

- Bilirubin is not stable in urine, especially when exposed to light.
- pH can affect urobilinogen levels. Alkaline urine indicates higher levels; acidic urine may show lower levels.
- ✚ Phenazopyridine colors the urine orange. This may give the false impression that the patient has jaundice.
- ✚ Cholestatic drugs may *decrease* urobilinogen levels.
- ✚ Antibiotics reduce intestinal flora, which in turn *decreases* urobilinogen levels.

### Crystals

- Radiographic contrast media may cause urinary crystals.

### WBCs

- Vaginal discharge may contaminate the urine specimen and factitiously cause WBCs in the urine.

### RBCs

- Strenuous physical exercise may cause RBC casts.
- Traumatic urethral catheterization may cause RBCs.
- Overaggressive anticoagulant therapy or bleeding disorders tend to cause RBCs in the urine without concomitant disease.

## Procedure and patient care

### Before

- PT Explain the procedure to the patient.

### During

- Collect a fresh urine specimen in a urine container.
- If the urine specimen contains vaginal discharge or bleeding, a clean-catch or midstream specimen will be needed. This requires meticulous cleaning of the urinary meatus with an antiseptic preparation to reduce contamination of the specimen by external organisms. The cleansing agent then must be completely removed or it will contaminate the specimen. For the *midstream collection*, use the following procedural steps:
  1. Have the patient begin urinating in a bedpan, urinal, or toilet and then stop urinating (this washes the urine out of the distal urethra).
  2. Correctly position a sterile urine container into which the patient voids 3 to 4 oz of urine.
  3. Cap the container.
  4. Allow the patient to finish voiding.

## 940 urinalysis

- For ketones, this test can be performed immediately after collection with a dipstick.
- With a Ketostix, dip the reagent into the urine specimen and remove it. Read the strip in 15 seconds by comparing it with the color chart.
- For urine specific gravity, a first-voided specimen is the best.
- For protein, the first-voided specimen is the best; however, occasionally a 24-hour urine collection is preferred.
- See inside front cover for Routine Urine Testing.

### After

- Transport the urine specimen to the laboratory promptly.
- If the specimen cannot be processed immediately, refrigerate it. If the urine cannot be tested within 2 hours of collection, a preservative should be used.
- If a 24-hour urine collection is requested, the specimen should be refrigerated, or a preservative should be used.
- Casts will break up as urine is allowed to sit. Urine examinations for casts should be performed on fresh specimens.

### Abnormal findings

#### Appearance and color

Bacteria	Dehydration
Pus	Overhydration
Red blood cells	Diabetes insipidus
Certain foods (e.g., beets, carrots)	Fever
Drug therapy (see Table 37)	Excessive sweating
Pathologic conditions (e.g., bleeding from the kidney)	Jaundice

#### Odor

Infection	Maple syrup urine disease
Ketonuria	
UTI	Phenylketonuria
Rectal fistula	Hepatic failure



## pH

## ▲ Increased levels

Respiratory alkalosis  
 Metabolic alkalosis  
 Urea-splitting bacteria  
 Vegetarian diet  
 Renal failure with inability  
 to form ammonia  
 Gastric suction  
 Vomiting  
 Diuretic therapy  
 Renal tubular acidosis  
 UTI

## ▼ Decreased levels

Metabolic acidosis  
 Diabetes mellitus  
 Diarrhea  
 Starvation  
 Respiratory acidosis  
 Sleep  
 Pyrexia

## Protein

## ▲ Increased levels

Nephrotic syndrome  
 Diabetes mellitus  
 Multiple myeloma  
 Preeclampsia  
 Glomerulonephritis  
 Congestive heart failure  
 Malignant hypertension  
 Polycystic kidney disease  
 Diabetic glomerulosclerosis

Amyloidosis  
 Systemic lupus  
 erythematosus  
 Goodpasture syndrome  
 Renal vein thrombosis  
 Heavy-metal poisoning  
 Galactosemia  
 Bacterial pyelonephritis  
 Nephrotoxic drug therapy  
 Bladder tumor

## Glucose

## ▲ Increased levels

Diabetes mellitus  
 Pregnancy  
 Renal glycosuria  
 Hereditary defects in metabolism of other reducing  
 substances (e.g., galactose, fructose, pentose)  
 Nephrotoxic chemicals (e.g., carbon monoxide, mercury,  
 lead)

## 942 urinalysis

### Specific gravity

#### ▲ Increased levels

- Dehydration
- Pituitary tumor or trauma that causes syndrome of inappropriate antidiuretic hormone (SIADH)
- Decrease in renal blood flow (as in heart failure, renal artery stenosis, or hypotension)
- Glycosuria and proteinuria
- Water restriction
- Fever
- Excessive sweating
- Vomiting
- Diarrhea
- X-ray contrast dye

#### ▼ Decreased levels

- Overhydration
- Diabetes insipidus
- Renal failure
- Diuresis
- Hypothermia
- Glomerulonephritis
- Pyelonephritis

### Leukocyte esterase

Possible UTI

### Nitrites

Possible UTI

### Ketones

- |   |                       |
|---|-----------------------|
| Uncontrolled diabetes mellitus            | After anesthesia      |
| Starvation                                | Prolonged vomiting    |
| Excessive aspirin ingestion               | Anorexia nervosa      |
| Ketoacidosis of alcoholism                | Fasting               |
| Febrile illnesses in infants and children | High-protein diets    |
| Weight reduction diets                    | Isopropanol ingestion |
|   | Dehydration           |

**Crystals**

Renal stone formation  
Drug therapy  
UTI

**Granular casts**

Acute tubular necrosis	Chronic lead poisoning
UTI	Reaction after exercise
Glomerulonephritis	Stress
Pyelonephritis	Renal transplant rejection
Nephrosclerosis	

**Fatty casts**

Nephrotic syndrome	Glomerulonephritis
Diabetic nephropathy	Chronic renal disease

**Epithelial casts**

Glomerulonephritis	Ethylene glycol intoxication
Eclampsia	Acute renal allograft rejection
Heavy-metal poisoning	

**Waxy casts**

Chronic renal disease	Glomerulonephritis
Chronic renal failure	Renal transplant rejection
Diabetic nephropathy	Nephrotic syndrome
Malignant hypertension	

**Hyaline casts**

Proteinuria	Glomerulonephritis
Fever	Pyelonephritis
Strenuous exercise	Congestive heart failure
Stress	Chronic renal failure

## 944 urinalysis

### Red blood cells and casts

#### ▲ Increased RBC levels

- Glomerulonephritis
- Interstitial nephritis
- Acute tubular necrosis
- Pyelonephritis
- Renal trauma
- Renal tumor
- Renal stones
- Cystitis
- Prostatitis
- Traumatic bladder catheterization

#### ▲ Increased RBC cast levels

- Glomerulonephritis
- Subacute bacterial endocarditis
- Renal infarct
- Goodpasture syndrome
- Vasculitis
- Sickle cell anemia
- Malignant hypertension
- Systemic lupus erythematosus

### White blood cells and casts

#### ▲ Increased WBC levels

- Bacterial infection in the urinary tract

#### ▲ Increased WBC cast levels

- Acute pyelonephritis
- Glomerulonephritis
- Lupus nephritis

notes

**urinary stone analysis** (Renal calculus analysis)**Type of test** Urine**Normal findings**

All urinary stones are pathologic.

**Test explanation and related physiology**

Urinary stone analysis is performed to identify the chemicals that make up a kidney stone and to treat any underlying disease that may have caused the stone formation. This information is also used to determine the most effective methods to diminish the chance of another stone.

Approximately 80% of stones are composed of calcium oxalate (CaOx) and calcium phosphate (CaP), 10% of struvite (magnesium ammonium phosphate produced during infection with bacteria that possess the enzyme urease), and 9% of uric acid (UA), and the remaining 1% are composed of cystine or ammonium acid urate or are diagnosed as drug-related stones.

Analysis is done on a kidney stone that has been passed in the urine or removed from the urinary tract during surgery. Analysis includes an evaluation of the size, shape, color, and weight of the stone.

**Interfering factors**

- Tape used to attach a stone to paper may affect the ability to accurately identify the composition of the stone.

**Procedure and patient care****Before**

**PT** Explain the procedure to the patient.

- Ensure pain relief if the patient is having ureteral colic.
- Obtain a history of any previous urinary stones.

**PT** Explain that there are no dietary restrictions for this test.

**PT** Provide and explain the use of the urine strainer.

**During**

**PT** Instruct the patient to urinate into the strainer provided.

**PT** Tell the patient to transfer any particulate matter to a container for laboratory analysis.

## 946 urinary stone analysis

### After

- Transport the specimen to the laboratory promptly.

### Abnormal findings

Urinary stones

notes

## urine culture and sensitivity (C&S)

**Type of test** Urine; microscopic

### Normal findings

Negative: < 10,000 bacteria per milliliter of urine

Positive: > 100,000 bacteria per milliliter of urine

### Test explanation and related physiology

Urine cultures and sensitivities are obtained to determine the presence of pathogenic bacteria in patients with suspected UTIs. An important part of any routine culture is to assess the *sensitivity* of any bacteria that are growing in the urine to various antibiotics. The physician can then more appropriately recommend the correct antibiotic therapy.

Most often, UTIs are limited to the bladder. However, the kidneys, ureters, bladder, or urethra can be the source of infection. All cultures should be performed before antibiotic therapy is initiated; otherwise the antibiotic may interrupt the growth of the organism in the laboratory. Most organisms require approximately 24 hours to grow in the laboratory, and a preliminary report can be given at that time. Usually 48 to 72 hours is required for growth and identification of the organism. With DNA sequencing through PCR or nanopore technology, the infecting bacteria can be identified and antibiotic sensitivity can be available within 4 hours of testing, allowing effective antibiotic therapy to be instigated very soon after testing.

To save money, urine cultures are usually done only if the urinalysis suggests a possible infection (e.g., increased number of WBCs, bacteria, high pH, positive leukocyte esterase).

### Interfering factors

- Contamination of the urine with stool, vaginal secretions, hands, or clothing will cause false-positive results.
- ✚ Drugs that may affect test results include antibiotics.

### Procedure and patient care

#### Before

- PT** Explain to the patient the procedure for obtaining a clean-catch (midstream) urine collection.
- Hold antibiotics until after specimen collection.
  - Provide the patient with the necessary supplies.

### During

- Note that a *clean-catch* or *midstream urine* collection is required for C&S testing. This requires meticulous cleansing of the urinary meatus with an antiseptic preparation to reduce contamination of the specimen by external organisms. Then the cleansing agent must be completely removed, or it will contaminate the urine specimen. The midstream collection is obtained by
  1. Having the patient begin to urinate in a bedpan, urinal, or toilet and then stop urinating (this washes the urine out of the distal urethra)
  2. Correctly positioning a sterile urine container, into which the patient voids 3 to 4 oz of urine
  3. Capping the container
  4. Allowing the patient to finish voiding
- Note that *urinary catheterization* may be needed for patients unable to void.
- For patients with an *indwelling urinary catheter*, obtain a specimen by aseptically inserting a syringe into the sampling port on the catheter. (Usually the catheter tubing distal to the puncture site needs to be clamped for 15-30 minutes before the aspiration of urine to allow urine to fill the tubing. After the specimen is withdrawn, the clamp is removed.)
- Note that *suprapubic aspiration* of urine is a safe method of obtaining urine in neonates and infants. The abdomen is prepared with an antiseptic, and a 25-gauge needle is inserted into the suprapubic area 1 inch above the symphysis pubis. Urine is aspirated into the syringe and then transferred to a sterile urine container.
- Collect specimens from infants and young children in a disposable pouch called a *U bag*. This bag has an adhesive backing around the opening to attach to the child.
- Note that for patients with a *urinary diversion* (e.g., an ileal conduit), catheterization should be done through the stoma.
- Urine should *not* be collected from an ostomy pouch.

### After

- Transport the specimen to the laboratory immediately (at least within 30 minutes).
- Notify the physician of any positive results so that appropriate antibiotic therapy can be initiated.

### Abnormal findings

UTI

notes



**urodynamic studies** (Uroflowmetry, Urine flow studies, Cystometry, Cystometrogram [CMG], Urethral pressure profile [UPP], Urethral pressure measurements)

**Type of test** Manometric

### Normal findings

Normal sensations of fullness and temperature

Normal pressures and volumes

Maximal cystometric capacity

Male: 350-750 mL

Female: 250-550 mL

Intravesical pressure when bladder is empty: usually < 40 cm H<sub>2</sub>O

Detrusor pressure: < 10 cm H<sub>2</sub>O

Residual urine: < 30 mL

### Maximal urethral pressures in normal patients

Age (years)	Male (cm H <sub>2</sub> O)	Female (cm H <sub>2</sub> O)
< 25	37-126	55-103
25-44	35-113	31-115
45-64	40-123	40-100
> 64	35-105	35-75

### Test explanation and related physiology

These tests are used to measure pressures and urine volume of the bladder and urethra in order to identify patients who have bladder function problems. Urodynamic tests usually include urine flow studies, postvoid residual (PVR) urine measurement, and CMG. They are used in patients with bladder outlet obstruction, urinary incontinence, and questionable neurogenic bladder. They are also used to document progress with treatment.

Both the motor and sensory function of the bladder are evaluated. During CMG, water is used to assess the first sensation of filling, fullness, and urinary urge. Bladder compliance and the presence of uninhibited detrusor contractions (i.e., phasic contractions) can also be noted during this filling CMG. Abdominal leak point pressure (ALPP) is useful in determining the etiology of urinary stress incontinence. The neuromuscular function of the bladder is assessed by measuring the efficiency of the detrusor muscle, intravesical pressure and capacity, and the bladder's response to thermal stimulation. Because urodynamic

studies have a wide range of normal, urodynamic findings of significance must be associated with reproduction of the patient's symptoms.

Cystometry can determine whether a bladder function abnormality is caused by neurologic, infectious, or obstructive diseases. Cystometry is indicated to elucidate the causes of bladder outlet obstruction or frequency and urgency. Cystometry is also part of the evaluation for the following: incontinence, persistent residual urine, vesicoureteral reflux, severe nocturnal enuresis, motor and sensory disorders affecting the bladder, and the effect of certain drugs on bladder function.

Uroflowmetry is the simplest of the urodynamic techniques, being noninvasive and requiring uncomplicated and relatively inexpensive equipment. It measures the volume of urine expelled from the bladder per second. If the rate is reduced, outflow obstruction can be documented and measured. Nomograms of maximal flow versus voided volume may be used for accurate test result interpretation, taking into account the patient's gender and age. Urine flowmeters provide a permanent graphic recording.

UPP is often performed with urine flow studies. UPP is the fluid pressure that would hypothetically be required to force open the collapsed urethra and allow urine to flow. The urethral pressure varies from point to point within the urethra. Thus, a graph of urethral pressure against distance along the urethra is determined. It is used to document reduced urethral pressures in incontinent patients. It is also used to indicate the degree of compression applied to the urethra from an abnormally enlarged prostate (which will increase UPP value).

Urodynamic studies are performed by a urologist in approximately 45 minutes and are often performed at the same time as cystoscopy. The only discomfort is that associated with the urethral catheterization.

### **Contraindications**

- UTIs because of the possibility of false results and the potential for the spread of infection.

### **Procedure and patient care**

#### **Before**

- PT** Explain the purpose and the procedure to the patient. Instruct the patient to arrive with a full bladder.
- Obtain informed consent.

- PT** Tell the patient that no fluid or food restrictions are needed.
- Assure the patient that he or she will be draped to prevent unnecessary exposure.
  - Assess the patient for signs and symptoms of UTI.
- PT** Instruct the patient not to strain while voiding because the results can be skewed.

### During

- Note the following procedural steps:
  1. Cystometry, usually performed in a urologist's office or a special procedure room, begins with the patient being asked to void.
  2. The amount of time required to initiate voiding and the size, force, and continuity of the urinary stream are recorded. The amount of urine, the time of voiding, and the presence of any straining, hesitancy, or terminal urine dribbling are also recorded.
  3. Next the bladder is tested for post void residual (see pelvic ultrasound, p. 685).
  4. A standing cough stress test is performed.
  5. The patient is placed in a lithotomy or supine position.
  6. A dual lumen catheter is inserted through the urethra and into the bladder.
  7. Residual urine volume is measured and recorded.
  8. Thermal sensation is evaluated by the instillation of approximately 30 mL of room-temperature saline solution into the bladder followed by an equal amount of warm water. The patient reports any sensations.
  9. This fluid is withdrawn from the bladder.
  10. The urethral catheter is connected to a cystometer to monitor bladder pressure. Sterile water, normal saline solution, or carbon dioxide gas is slowly introduced into the bladder at a controlled rate, usually with the patient in a sitting position. While the bladder is slowly filled, pressures are simultaneously recorded. This is called a cystometrogram.
  11. Patients are asked to indicate the first urge to void and then when they have the feeling that they must void. The bladder is full at this point.
  12. The pressures and volumes are plotted on a graph.
  13. ALPP to investigate for stress urinary incontinence are obtained by asking the patient to perform the Valsalva maneuver in gradients (i.e., mild, moderate, strong) followed by cough (i.e., mild, moderate, strong).

## 952 urodynamic studies

14. The patient is asked to void around the catheter, and the maximal intravesical voiding pressure is recorded.
  15. The bladder is drained for any residual fluid or gas. If no additional studies are to be done, the urethral catheter is removed.
  16. For urethral pressures, fluid or gas is instilled through the catheter, which is withdrawn while pressures along the urethral wall are obtained
- Note that pelvic floor sphincter electromyography (see p. 353) can be performed to evaluate the urethral sphincter in cases of incontinence.
- PT** Throughout the study, ask the patient to report any sensations, such as pain, flushing, sweating, nausea, bladder filling, and an urgency to void.
- Note that certain drugs may be administered during the cystometric examination to distinguish between underactivity of the bladder because of muscle failure and underactivity associated with denervation. Cholinergic drugs (e.g., bethanechol [Urecholine]) may be given to enhance the tone of a flaccid bladder. Anticholinergic drugs (e.g., atropine) may be given to promote relaxation of a hyperactive bladder. If these drugs are to be given, the catheter is left in place. After these drugs are given, the examination is repeated 20 to 30 minutes later using the first test as a control value. The information obtained with the drugs assists in deciding whether drugs will be effective treatment.

### After

- Observe the patient for any manifestations of infection (e.g., elevated temperature, chills, or dysuria).
- Examine the urine for hematuria. Notify the physician if the hematuria persists after several voidings.
- Provide a warm sitz bath or tub bath for the patient's comfort if desired.

### Abnormal findings

Neurogenic bladder

Bladder obstruction, infection, or hypertonicity

Diminished bladder capacity

Prostatic obstruction secondary to benign prostatic hypertrophy or cancer

Urinary incontinence

### notes

## uroporphyrinogen-1-synthase

**Type of test** Blood

### Normal findings

1.27-2 mU/g of hemoglobin or 81.9-129.6 units/mol Hgb (SI units)

### Test explanation and related physiology

Porphyria is a group of genetic disorders characterized by an accumulation of porphyrin products, usually in the liver. This group of disorders results from enzymatic deficiencies in synthesis of heme (a part of hemoglobin). Acute intermittent porphyria is the most common form of the liver porphyrias; it is caused by a deficiency in uroporphyrinogen-1-synthase (also called porphobilinogen deaminase). This enzyme is necessary for erythroid cells to make heme.

This enzyme is significantly reduced during the acute and latent phases of this disorder. It is important to identify this deficiency to prevent acute bouts of porphyria. See porphyrins and porphobilinogens on p. 717.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: purple
- Because this test is based on the hemoglobin measurement, concurrently obtain a hemoglobin level on the patient.
- Indicate on the laboratory slip if the patient is having symptoms of acute porphyria.

### Abnormal findings

- ▼ Decreased levels  
Acute intermittent porphyria

notes

**vanillylmandelic acid** (VMA), homovanillic acid (HVA), and catecholamines (Epinephrine, Norepinephrine, Metanephrine, Normetanephrine, Dopamine)

**Type of test** Urine (24-hour)

### Normal findings

#### VMA

Adult/elderly: < 6.8 mg/24 hr or < 35  $\mu$ mol/24 hr (SI units)

Adolescent: 1-5 mg/24 hr

Child: 1-3 mg/24 hr

Infant: < 2 mg/24 hr

Newborn: < 1 mg/24 hr

#### HVA

$\geq$ 15 years (adults): not applicable

10-14 years: < 12 mg/g creatinine

5-9 years: < 9 mg/g creatinine

2-4 years: < 13.5 mg/g creatinine

1 year: < 23 mg/g creatinine

< 1 year: 35 mg/g creatinine

#### Catecholamines

Free catecholamines

< 100 mcg/24 hr or < 590 nmol/day (SI units)

Epinephrine

Adult/elderly: < 20 mcg/24 hr or < 109 nmol/day (SI units)

Child

0-1 year: 0-2.5 mcg/24 hr

1-2 years: 0-3.5 mcg/24 hr

2-4 years: 0-6 mcg/24 hr

4-7 years: 0.2-10 mcg/24 hr

7-10 years: 0.5-14 mcg/24 hr

Norepinephrine

Adult/elderly: < 100 mcg/24 hr or < 590 nmol/day (SI units)

Child

0-1 year: 0-10 mcg/24 hr

1-2 years: 0-17 mcg/24 hr

2-4 years: 4-29 mcg/24 hr

4-7 years: 8-45 mcg/24 hr

7-10 years: 13-65 mcg/24 hr

Dopamine

Adult/elderly: 65-400 mcg/24 hr

## Child

0-1 year: 0-85 mcg/24 hr

1-2 years: 10-140 mcg/24 hr

2-4 years: 40-260 mcg/24 hr

&gt; 4 years: 65-400 mcg/24 hr

## Metanephrine

< 1.3 mg/24 hr or < 7  $\mu$ mol/day (SI units)

## Normetanephrine

15-80 mcg/24 hr or 89-473 nmol/day (SI units)

**Test explanation and related physiology**

This 24-hour urine test is a screening test for the diagnosis of catecholamine-producing tumors such as neuroblastoma, pheochromocytoma, and other rare adrenal and neural crest tumors. Likewise, neural crest tumors such as neuroblastoma can also hypersecrete catecholamines. Dopamine is the precursor of epinephrine and norepinephrine. HVA is a metabolite of dopamine. Metanephrine and normetanephrine are catabolic products of epinephrine and norepinephrine, respectively. VMA (3-methoxy-4-hydroxymandelic acid) is the product of catabolism of both metanephrine and normetanephrine. In pheochromocytoma, one or all of these substances will be present in excessive quantities in a 24-hour urine collection. These hormones may be measured singularly in the urine, but the collective metabolic end products, HVA and VMA, are more easily detected because their concentrations are much higher than any one catecholamine component.

VMA and HVA are primarily used as a screening test for neural crest tumors. These urinary tests can also be used to monitor tumor activity. (See also metanephrine, plasma, p. 621.)

**Interfering factors**

- Increased levels of VMA may be caused by certain foods (e.g., tea, coffee, cocoa, vanilla, chocolate).
- Vigorous exercise, stress, and starvation may cause increased VMA levels.
- Falsely decreased levels of VMA may be caused by uremia, alkaline urine, and radiographic iodine contrast agents.
- ✚ Drugs that may cause *increased* VMA levels include caffeine, epinephrine, levodopa, lithium, and nitroglycerin.
- ✚ Drugs that may cause *decreased* VMA levels include clonidine, disulfiram, guanethidine, imipramine, monoamine oxidase inhibitors, phenothiazines, and reserpine.

## 956 vanillylmandelic acid

- ❖ Drugs that may cause *increased* catecholamine levels include alcohol (ethyl), aminophylline, caffeine, chloral hydrate, clonidine (chronic therapy), contrast media (iodine containing), disulfiram, epinephrine, erythromycin, insulin, methenamine, methyldopa, nicotinic acid (large doses), nitroglycerin, quinidine, riboflavin, and tetracyclines.
- ❖ Drugs that may cause *decreased* catecholamine levels include guanethidine, reserpine, and salicylates.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient.
- PT** Explain the dietary restrictions and the 24-hour urine collection procedure to the patient.
- PT** For 2 or 3 days before the 24-hour collection for VMA and throughout the collection, place the patient on a VMA-restricted diet. Generally, instruct the patient to avoid coffee, tea, bananas, chocolate, cocoa, licorice, citrus fruit, all foods and fluids containing vanilla, and aspirin. Obtain specific restrictions from the laboratory.
- PT** Inform the patient of the need to avoid taking antihypertensive medications, and sometimes all medications, during this period and possibly even longer.

#### During

- Collect the 24-hour urine specimen using a preservative.
- PT** Instruct the patient to begin the 24-hour urine collection after voiding. See inside front cover for Routine Urine Testing.
- Identify and minimize factors contributing to patient stress and anxiety. Excessive physical exercise and emotion may alter catecholamine test results by causing an increased secretion of epinephrine and norepinephrine.

#### After

- Send the specimen to the laboratory as soon as the test is completed.
- Allow the patient to have foods and drugs that were restricted in preparation for the test.



## Abnormal findings

### ▲ Increased levels

- Pheochromocytomas
- Neuroblastomas
- Ganglioneuromas
- Ganglioblastomas
- Severe stress
- Strenuous exercise
- Acute anxiety

### notes

**varicella virus testing** (Antibody testing for Varicella zoster virus, Varicella zoster virus [VZV], Herpes virus-3)**Type of test** Microscopic**Normal findings**IgM  $\leq 0.90$  ISR (immune status ratio)

IgG

Vaccinated: positive ( $\geq 1.1$  AI [antibody index])Unvaccinated: negative ( $\leq 0.8$  AI)

Culture: no growth

**Test explanation and related physiology**

Varicella is an acute infectious disease caused by varicella zoster virus (VZV). The recurrent infection, herpes zoster, is also known as shingles. VZV is a DNA virus and is a member of the herpesvirus group that includes herpes simplex types 1 and 2; see p. 498). Like other herpesviruses, VZV has the capacity to persist in the body after the primary (first) infection as a recurrent and latent infection. VZV persists in sensory nerve ganglia. Primary infection with VZV results in chickenpox. The recurrent or latent infections cause herpes zoster (shingles).

Testing for VZV or the antibodies produced in response to VZV is not routinely used to diagnose active cases of chickenpox and shingles, which are detected by the person's signs and symptoms. It may be performed in pregnant women, newborns, in people before organ transplantation, and in those with HIV/AIDS. Testing may also be used to determine if someone has been previously exposed to VZV either through past infection or vaccination and has developed immunity to the disease. It can distinguish between an active or prior infection. It can determine whether someone with severe or atypical symptoms has an active VZV infection or has another condition with similar symptoms.

Recovery from primary varicella infection usually results in lifetime immunity. Herpes zoster, or shingles, occurs when latent VZV reactivates and causes recurrent disease. Postherpetic neuralgia (pain in the area of the zoster) may persist after the zoster has resolved. This can be reduced by aggressive and early antiviral treatment.

Testing is also recommended to confirm varicella outbreaks or establish varicella as a cause of death. A variety of serologic tests for varicella antibody are available commercially and can be used to assess disease-induced immunity. The most frequent source of

VZV isolation is vesicular fluid. Serum is used for serologic testing and some molecular assays.

Live-attenuated varicella virus is now part of the measles–mumps–rubella and varicella vaccine (e.g., ProQuad) for use in children. A herpes zoster vaccine is available for persons 50 years of age and older to minimize the risk of herpes zoster.

## Procedure and patient care

### Before

**PT** Explain the procedure to the patient.

**PT** Tell the patient that no preparation is required.

### During

- Specimens are best collected by unroofing a vesicle, preferably a fresh fluid-filled vesicle, and then rubbing the base of a skin lesion with a polyester swab.
- For serologic testing, collect whole venous peripheral blood in serum separator (gold, marble or grey) vacutainer.
- *For the blood spot method:*
  - Prick the subject's finger using a lancet.
  - Collect a sufficient quantity of blood onto both defined areas on the filter strip so that the spot expands to the circular border. Permit the specimen to air dry completely.

### After

- Check for bleeding.

## Abnormal findings

Acute varicella infection

Herpes (varicella) zoster

Varicella immunity

### notes

**vascular ultrasound studies** (Venous Doppler ultrasound, Venous duplex scan)

**Type of test** Ultrasound

**Normal findings****Venous**

A normal Doppler venous signal with spontaneous respiration  
Normal venous system without evidence of occlusion

**Arterial**

Normal arterial Doppler signal with systolic and diastolic components

No reduction in blood pressure in excess of 20 mm Hg compared with the normal extremity

A normal ankle-to-brachial arterial blood pressure index of 0.85 or greater

No evidence of arterial occlusion

**Test explanation and related physiology**

Vascular ultrasound studies are used to identify occlusion or thrombosis of the veins or arteries of an extremity. Patency is demonstrated with Doppler ultrasound by detecting moving red blood cells (RBCs) in the vessel. The patency of the venous system can also be identified by evaluating the degree of venous reflux (backward blood flow in the veins of the lower extremities in patients with venous valvular insufficiency).

Vascular duplex scanning combines the benefits of Doppler with B-mode scanning. With this technique, one is able to directly visualize areas of vascular narrowing or occlusion. The degree of occlusion is measured as percentage of the entire lumen that is occluded. Venous thrombosis is suspected when the vein is not easily compressible by the ultrasound probe.

*Color Doppler ultrasound (CDU)* can be added to arterial duplex scanning. This allows visualization of stenotic areas based on velocity or direction of blood flow in a particular area of the artery.

Duplex scanning is routinely used to identify venous thrombosis in patients suspected of having deep vein thrombosis (DVT) in an extremity. It is more rapidly performed and interpreted than venography (see p. 963) or MRV (see p. 601). In general, venous duplex scanning is less accurate than venography in identifying DVT in the calf or in the iliac veins.

With a single-mode transducer, blood flow can be heard audibly and is augmented by an audio speaker as a swishing noise. With single-mode arterial Doppler studies, peripheral arteriosclerotic occlusive disease of the extremities can be easily located. By slowly deflating blood pressure cuffs placed on the calf and ankle, systolic pressure in the arteries of the extremities can be accurately measured by detecting the first evidence of blood flow with the Doppler transducer. The extremely sensitive Doppler ultrasound detector can recognize the swishing sound of even the most minimal blood flow. Normally systolic blood pressure is slightly higher in the arteries of the arms than in the legs. If the difference in blood pressure exceeds 20 mm Hg, occlusive disease is believed to exist immediately proximal to the area tested. Lower extremity arterial bypass graft patency can also be assessed with Doppler ultrasound.

### Interfering factors

- Venous or arterial occlusive disease proximal to the site of testing
- Cigarette smoking because nicotine can cause constriction of the peripheral arteries and alter the results

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

**PT** Inform the patient that this is a painless procedure.

- Remove all clothing from the extremity to be examined.

**PT** Instruct the patient to abstain from cigarette smoking for at least 30 minutes before the test.

#### During

- Note the following procedural steps:

##### *Venous studies*

1. A conductive gel is applied to the skin overlying the venous system of the extremity in multiple areas.
2. Usually, for the lower extremity, the deep venous system is identified in the ankle, calf, thigh, and groin.
3. The characteristic “swishing” sound of the blood flow indicates a patent venous system.
4. Usually both the superficial and deep venous systems are evaluated.

##### *Arterial studies*

1. These are performed with the use of blood pressure cuffs, which are placed around the thigh, calf, and ankle.

## 962 vascular ultrasound studies

2. A conductive paste is applied to the skin overlying the artery distal to the cuffs.
  3. The proximal cuff is inflated to a level higher than systolic blood pressure in the normal extremity.
  4. The Doppler ultrasound transducer is placed immediately distal to the inflated cuff.
  5. The pressure in the cuff is slowly released.
  6. The highest pressure at which blood flow is detected by the characteristic “swishing” Doppler signal is recorded as the blood pressure of that artery.
  7. The test is repeated at each successive level.
  8. When the ankle pressure is divided by the arm (brachial artery) pressure, this is known as the ankle brachial index (ABI). If the ABI is less than 0.85, significant arterial occlusive disease exists within the extremity.
- Note that these studies are usually performed in the vascular laboratory or radiology department and take approximately 30 minutes.

### After

- Encourage the patient to verbalize his or her concerns.
- Remove the transducer gel from the extremity.

**PT** Inform the patient that the physician must interpret the studies and that results will be available in a few hours.

### Abnormal findings

Venous occlusion secondary to thrombosis or thrombophlebitis

Venous varicosities

Small or large vessel arterial occlusive disease

Spastic arterial disease (e.g., Raynaud phenomenon)

Small vessel arterial occlusive disease (as in diabetes)

Embolic arterial occlusion

Arterial aneurysm

### notes

## venography of lower extremities (Phlebography, Venography)

**Type of test** X-ray with contrast dye

### Normal findings

No evidence of venous thrombosis or obstruction

### Test explanation and related physiology

Venography is an x-ray study designed to identify and locate thrombi within the venous system (most commonly in the extremities). During this study, dye is injected into the venous system of the affected extremity. X-ray images are then taken at timed intervals to visualize the venous system. Obstruction to the flow of dye or a filling defect within the dye-filled vein indicates that thrombosis exists. A positive study accurately confirms the diagnosis of venous thrombosis; however, a normal study, although not as accurate, does make the diagnosis of venous thrombosis very unlikely. Venography is also used to identify venous stenosis caused by external obstruction or indwelling catheter-induced thrombosis.

Often both extremities are studied, even though only one leg is suspected of containing deep vein thrombosis. The normal extremity is used for comparison with the involved extremity. Venography is more accurate than venous Doppler (see p. 960) for thrombosis below the knee or in the femoral veins.

### Contraindications

- Patients with severe edema of the legs
- Patients who are uncooperative
- Patients who are allergic to iodinated dye
- Patients with renal failure because the dye is nephrotoxic

### Potential complications

- Allergic reaction to iodinated dyes
- Renal failure, especially in elderly people who are chronically dehydrated or may have a mild degree of renal failure
- Subcutaneous dye infiltration, causing cellulitis and pain
- Venous thrombophlebitis caused by the dye
- Bacteremia caused by a break in a sterile technique
- Venous embolism caused by dislodgment of a deep vein clot induced by the dye injection
- Lactic acidosis may occur in patients who are taking metformin and receiving iodine contrast. The metformin should be held the day of the test to prevent this complication.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Obtain informed consent for this procedure.
  - Assess the patient for allergies to iodinated dyes and shellfish.
  - If needed, provide appropriate pain medication.
- PT** Ensure that the patient is appropriately hydrated before testing. Injection of the iodinated contrast may cause renal failure, especially in elderly patients.

### During

- Note the following procedural steps:
  1. The patient is taken to the radiology department and placed in a supine position on the x-ray table.
  2. Catheterization of a superficial vein on the foot is performed. This may require a surgical cutdown.
  3. An iodinated, radiopaque dye is injected into the vein.
  4. X-ray images are taken to follow the course of the dye up the leg.
  5. A tourniquet is frequently placed on the leg to prevent filling of the superficial saphenous vein. As a result, all of the dye goes toward filling the deep venous system, which contains the most clinically significant thrombosis that can embolize.
- Note that a radiologist performs this study in approximately 30 to 90 minutes.
- PT** Tell the patient that the venous catheterization is only as uncomfortable as a cutaneous heel stick.
- PT** The dye may cause the patient to feel a warm flush. Occasionally, mild degrees of nausea, vomiting, or skin itching also may occur.

### After

- Continue appropriate hydration of the patient.
- Observe the puncture site for infection, cellulitis, or bleeding.
- Assess the patient's vital signs for signs of bacteremia.

## Abnormal findings

Obstructed venous systems

Acute deep vein thrombosis

notes



## viral cultures

**Type of test** Blood, urine, stool, throat, skin lesions

### Normal findings

No virus isolated

### Test explanation and related physiology

Viral infections are the most common infections affecting children and adults. Viruses are subdivided by the nucleic acid material they contain (RNA or DNA). Infections from viruses are often indistinguishable from bacterial infections. Definitive diagnosis of viral disease is made by culture of the virus (discussed here). Other methods used to identify viral disease include:

- Serologic methods of identifying antibodies to a specific virus
- Serologic methods of identifying antigen parts of a virus
- Direct detection by electron microscopy
- Detection by nucleic acid probes/viral load

The ability to isolate a virus in culture depends on many aspects of the culture process. The first is determining the correct specimen for culture. That depends on the organ involved and the type of virus suspected ([Table 38](#)). Timing is important. Viral load is always the greatest in the early stages of the disease. This can now be measured by direct quantification of viral DNA or RNA. Cultures obtained in the first few days after symptoms appear offer the best chance of identifying the infective culture. Viral cultures take 1 to 2 days to be reported.

**TABLE 38** Specimen culture for common viruses and diseases

Common viruses	Specimens	Diseases
Adenovirus, influenza, respiratory syncytial, rhinovirus	Throat culture, bronchoscopic aspiration, nasopharyngeal swab	Influenza, pneumonia, pharyngitis
Rubella, rubeola, coxsackie, varicella	Throat culture, skin vesicle	Skin rash, zoster
Arbovirus, enterovirus, herpes simplex, cytomegalovirus	Throat culture, cerebrospinal fluid, blood	Meningitis, encephalitis
Influenza A	Throat culture	Flu syndrome

## Interfering factors

- An inadequate specimen, the timing, or the choice of culture medium will cause false-negative test results.
- The use of a cotton swab or wooden applicator for specimen collection may destroy the virus.

## Procedure and patient care

### Before

- PT** Explain the method of collection of the specimen to the patient.
- Obtain a history regarding the timing of symptoms.
- Accurately record the source of the specimen.

### During

- Use a closed specimen system to obtain and transport the specimen to the laboratory.
- Transport the specimen immediately to the laboratory. Viruses in specimens quickly lose their vitality.
- Place samples on ice if delivery to the laboratory is not immediate.
- Small-volume specimens such as tissue aspirates are often best transported in viral transport medium. If bacterial cultures are to be performed, collect as a separate specimen.
- If blood is the specimen, obtain 5 mL to 7 mL of blood in a lavender- or blue-top tube.

### After

- PT** Explain that the patient may still be infectious and should minimize exposure to others.

## Abnormal findings

Viral infectious disease (see [Table 38](#))

### notes

## virus testing

**Type of test** Blood, miscellaneous

### Normal findings

Negative for viral antibody

### Test explanation and related physiology

Testing for a virus is indicated when a person with viral symptoms lives in or has traveled to an area harboring the virus. Testing is done in the clinical setting when a patient has severe symptoms contributing to significant morbidity. Testing is also performed for epidemiologic reasons to identify a viral outbreak and its extent. Finally, testing can indicate immunity after exposure to the virus or a vaccination.

Viral testing is performed by identifying antibodies to viral particles in the blood or in other body fluids. Virus can be cultured in special media and subsequently identified (see viral culture, p. 965). Finally viral RNA or DNA can be identified in body fluids (e.g., nasopharyngeal mucus or blood) heavily contaminated by the live virus. Epstein-Barr (see p. 368), hepatitis (see p. 493), respiratory syncytial, herpes, parainfluenza, HIV (see p. 506), Dengue fever, coxsackie, choriomeningitis, mumps, West Nile virus (WNV), arbovirus, equine, cytomegalovirus, rubella, and influenza A and B (including  $H_1N_1$ ) can all be detected by identifying antibodies to some part of the viral particle. Testing is now available to identify antibodies to many viruses in the same sample, indicating all past and present viral infections affecting a particular patient.

Front-line testing measures IgM or IgG antibodies to the virus that may or may not be specific to that particular virus. If the front-line testing is positive, confirmatory tests may be carried out. This testing may be particularly important for public health officials and researchers. Some viruses (e.g., WNV) can be transmitted through donated blood or blood components. For that reason, WNV testing kits for WNV antibodies are routinely performed on all donated blood.

### Interfering factors

- Other viral infections will cause elevations of serologic testing, especially when combined total immunoglobulin (Ig) M and IgG are tested.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red
- Obtain other samples for testing as directed. These may include nasopharyngeal swabs, sputum, or other body fluids.
- PT** Explain to the patient and family that, in most circumstances, testing is carried out at a referral laboratory.
- PT** Inform patients that they should observe isolation precautions until results are negative. Explain that results may not be available for 2 weeks.

### Abnormal findings

Viral infections (acute and chronic)

notes

**vitamin B<sub>12</sub>** (Cyanocobalamin and methylmalonic acid [MMA])**Type of test** Blood; urine**Normal findings**Vitamin B<sub>12</sub>: 160-950 pg/mL or 118-701 pmol/L (SI units)

MMA: &lt; 3.6 μmol/mmol creatinine

**Test explanation and related physiology**

Vitamin B<sub>12</sub> is necessary for conversion of the inactive form of folate to the active form. This function is most notable in the formation and function of RBCs. Vitamin B<sub>12</sub> deficiency, like folic acid deficiency, causes anemia. The RBCs formed in light of these deficiencies consist of large megaloblastic RBCs. These RBCs cannot conform to the size of small capillaries. Instead, they fracture and hemolyze. The shortened life span ultimately leads to anemia.

In the stomach, gastric acid detaches vitamin B<sub>12</sub> from its binding proteins. Intrinsic factor (IF), which is necessary for vitamin B<sub>12</sub> absorption in the terminal ileum, is made in the stomach mucosa. Without IF, vitamin B<sub>12</sub> cannot be absorbed. Deficiency of IF is the most common cause of vitamin B<sub>12</sub> deficiency (pernicious anemia). The next most common cause of vitamin B<sub>12</sub> deficiency is lack of gastric acid to separate the ingested vitamin B<sub>12</sub> from its binding proteins. This is common in patients who have had gastric surgery. A third cause of vitamin B<sub>12</sub> deficiency is malabsorption caused by diseases of the small terminal ileum. Vitamin B<sub>12</sub> deficiency is common in poorly nourished elderly people and in vegetarians.

Serum B<sub>12</sub> is a measurement of recent B<sub>12</sub> ingestion. More prolonged B<sub>12</sub> deficiency is better and more easily measured by *urinary methylmalonic acid (MMA)* measurement. Elevated serum MMA levels and urinary excretion of MMA are direct measures of tissue vitamin B<sub>12</sub> activity. The active form of B<sub>12</sub> is essential in the intracellular conversion of L-methylmalonyl coenzyme A (MMA CoA) to succinyl CoA. Without B<sub>12</sub>, MMA CoA metabolism is diverted to make large quantities of MMA. MMA is then excreted by the kidneys. MMA testing is the most sensitive test for B<sub>12</sub> deficiency.

**Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: verify with laboratory
- Blood tube commonly used: red

## 970 vitamin B<sub>12</sub>

**PT** Instruct the patient not to consume alcoholic beverages before the test. Check the time period with the physician or laboratory.

- Draw the specimen before starting vitamin B<sub>12</sub> therapy.

### Abnormal findings

#### ▲ Increased levels

Leukemia  
Polycythemia vera  
Severe liver dysfunction  
Myeloproliferative disease

#### ▼ Decreased levels

Pernicious anemia  
Malabsorption syndromes  
Inflammatory bowel disease  
Intestinal worm infestation  
Atrophic gastritis  
Zollinger-Ellison syndrome  
Large proximal gastrectomy  
Resection of terminal ileum  
Achlorhydria  
Pregnancy  
Vitamin C deficiency  
Folic acid deficiency

notes

## **vitamin D** (25-hydroxy vitamin D<sub>2</sub> and D<sub>3</sub>; 1,25-dihydroxyvitamin D [1,25(OH)<sub>2</sub>D])

**Type of test** Blood

### **Normal findings**

Total 25-hydroxy D (D<sub>2</sub> + D<sub>3</sub>): 25-80 ng/mL

1,25(OH)<sub>2</sub>D

Males: 18-64 pg/mL

Females: 18-78 pg/mL

### **Test explanation and related physiology**

Vitamin D levels are calculated to ensure that postmenopausal women have adequate vitamin D levels to absorb dietary calcium. Because of the increased number of research studies investigating the role of vitamin D in osteoporosis and cancer prevention, more and more patients are having this blood test.

Vitamin D is a fat-soluble vitamin. The two major forms of vitamin D are vitamin D<sub>2</sub> (or ergocalciferol) and vitamin D<sub>3</sub> (or cholecalciferol). Vitamin D<sub>2</sub> is provided by dietary sources. Because only fish is naturally rich in vitamin D, most of the vitamin D<sub>2</sub> intake in the industrialized world is from fortified products, including milk, soy milk, and breakfast cereals or supplements.

Vitamin D<sub>3</sub> is produced in skin exposed to sunlight, specifically ultraviolet B (UVB) radiation. In this scenario, 7-dehydrocholesterol reacts with UVB light at wavelengths between 270 nm and 300 nm to produce vitamin D<sub>3</sub>.

After vitamin D is produced in the skin or consumed in food, it is converted in the liver and kidney to form 1,25-dihydroxyvitamin D (1,25[OH]<sub>2</sub>D), the physiologically active form of vitamin D. After this conversion, the hormonally active form of vitamin D is released into the circulation. Vitamin D regulates the calcium and phosphorus levels in the blood by promoting their absorption from food in the intestines and by promoting reabsorption of calcium in the kidneys. This enables normal mineralization of bone needed for bone growth and bone remodeling. Vitamin D inhibits parathyroid hormone secretion from the parathyroid gland. Vitamin D promotes the immune system by increasing phagocytosis, antitumor activity, and other immunomodulatory functions.

Vitamin D deficiency can result from inadequate dietary intake, inadequate sunlight exposure, malabsorption syndromes, liver or kidney disorders, or by a number of metabolic hereditary

disorders. Deficiency results in impaired bone mineralization and leads to bone softening diseases (rickets in children and osteomalacia in adults). Vitamin D deficiency may also contribute to the development of osteoporosis.

Vitamin D levels can be measured in the blood. Usually 25-hydroxy vitamin D<sub>2</sub> and D<sub>3</sub> are measured and added to obtain the total 25-hydroxy vitamin D level. They are usually reported individually and as a total. Therapy is based on the measurement of total 25-hydroxy vitamin D levels. 1,25-dihydroxyvitamin D (the active metabolite of vitamin D) can be measured and is helpful in patients who have signs of vitamin D deficiency yet have normal levels of total vitamin D.

### **Interfering factors**

- ✚ Corticosteroid drugs can *decrease* vitamin D levels by reducing calcium absorption.
- ✚ The weight-loss drug orlistat and the cholesterol-lowering drug cholestyramine can *decrease* vitamin D levels by reducing the absorption of vitamin D and other fat-soluble vitamins.
- ✚ Barbiturates and phenytoin *decrease* vitamin D levels by increasing hepatic metabolism of vitamin D to inactive compounds.

### **Procedure and patient care**

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: red or green
- PT** If the patient has a vitamin D deficiency, educate him or her about dietary food sources and about the importance of sunlight.



**Abnormal findings****▲ Increased levels**

- Williams syndrome
- Excess dietary supplements

**▼ Decreased levels**

- Rickets
- Osteomalacia
- Osteoporosis
- Gastrointestinal malabsorption syndromes
- Renal failure
- Liver disease
- Familial hypophosphatemic rickets (X-linked hypophosphatemic rickets)
- Acute inflammatory disease
- Inadequate dietary intake
- Inadequate exposure to sunlight

notes

## white blood cell count and differential count (WBC and differential, Leukocyte count)

**Type of test** Blood

### Normal findings

#### Total WBCs

Adult/child > 2 years: 5000-10,000/mm<sup>3</sup> or 5-10 × 10<sup>9</sup>/L  
(SI units)

Child ≤ 2 years: 6200-17,000/mm<sup>3</sup>

Newborn: 9000-30,000/mm<sup>3</sup>

#### Differential count

	%	Absolute (per mm <sup>3</sup> )
Neutrophils	55-70	2500-8000
Lymphocytes	20-40	1000-4000
Monocytes	2-8	100-700
Eosinophils	1-4	50-500
Basophils	0.5-1	25-100

### Possible critical values

WBC count < 2500 or > 30,000/mm<sup>3</sup>

### Test explanation and related physiology

The WBC count has two components. The first is a count of the total number of WBCs (leukocytes) in 1 mm<sup>3</sup> of peripheral venous blood. The other component, the differential count, measures the percentage of each type of leukocyte present in the same specimen. An increase in the percentage of one type of leukocyte means a decrease in the percentage of another. Neutrophils and lymphocytes make up 75% to 90% of the total leukocytes. The total leukocyte count has a wide range of normal values, but many diseases may induce abnormal values.

An increased total WBC count (leukocytosis: WBC count > 10,000/mm<sup>3</sup>) usually indicates infection, inflammation, tissue necrosis, or leukemic neoplasia. Trauma or stress, either emotional or physical, may increase the WBC count. A decreased total WBC count (leukopenia: WBC count < 4000/mm<sup>3</sup>) occurs in many forms of bone marrow failure (e.g., after antineoplastic chemotherapy or radiation therapy, marrow infiltrative diseases, overwhelming infections, dietary deficiencies, and autoimmune diseases).

The major function of the WBCs is to fight infection and react against foreign bodies or tissues. Five types of WBCs may easily be identified on a routine blood smear (see p. 149, blood smear). These cells, in order of frequency, include neutrophils, lymphocytes, monocytes, eosinophils, and basophils.

WBCs are divided into granulocytes and nongranulocytes. Granulocytes include neutrophils, basophils, and eosinophils. Neutrophils have multilobed nuclei and are sometimes referred to as polymorphonuclear leukocytes (PMNs or *polys*). The normal ranges for absolute counts depend on age, sex, and ethnicity. For example, the normal range for absolute neutrophils for African American men is 1400 to 7000 cells/microliter.

*Neutrophils*, the most common granulocyte, are produced in 7 to 14 days, and exist in the circulation for only 6 hours. The primary function of the neutrophil is phagocytosis (killing and digestion of bacterial microorganisms). Acute bacterial infections and trauma stimulate neutrophil production, resulting in an increased WBC count. Often when neutrophil production is significantly stimulated, early immature forms of neutrophils enter the circulation. These immature forms are called *band* or *stab* cells. This occurrence, referred to as a *shift to the left* in WBC production, is indicative of an ongoing acute bacterial infection.

*Basophils* (also called mast cells), especially *eosinophils*, are involved in the allergic reaction. Parasitic infestations also are capable of stimulating the production of these cells. These cells are capable of phagocytosis of antigen–antibody complexes. As the allergic response diminishes, the eosinophil count decreases. Eosinophils and basophils do not respond to bacterial or viral infections.

Nongranulocytes (mononuclear cells) include lymphocytes, monocytes, and histiocytes. *Lymphocytes* are divided into two types: T-cells and B-cells. Whereas T-cells are primarily involved with cellular-type immune reactions, B-cells participate in humoral immunity (antibody production). The primary function of the lymphocytes is fighting chronic bacterial and acute viral infections. The differential count does not separate the T- and B-cells but rather counts the combination of the two.

*Monocytes* are phagocytic cells capable of fighting bacteria in a way very similar to that of neutrophils. However, monocytes can be produced more rapidly and can spend a longer time in the circulation than neutrophils.

The WBC count and differential counts are routinely measured as part of the complete blood count (see p. 269). Serial WBC counts and differential counts have both diagnostic and

## 976 white blood cell count and differential count

prognostic value. For example, a persistent increase in the WBC count (particularly the neutrophils) may indicate a worsening of an infectious process (e.g., appendicitis). A dramatic decrease in the WBC count below the normal range may indicate marrow failure. In patients receiving chemotherapy, a reduced WBC count may delay further chemotherapy.

The absolute count is calculated by multiplying the differential count (%) by the total WBC count. For example, the *absolute neutrophil count* (ANC) is helpful in determining the patient's real risk for infection. It is calculated by multiplying the WBC count by the percent of neutrophils and percent of bands; that is,

$$\text{ANC} = \text{WBC} \times (\% \text{ neutrophils} + \% \text{ bands})$$

If the ANC is less than 1000, the patient may need to be placed in protective isolation because he or she could be severely immunocompromised and at great risk for infection.

### Interfering factors

- Physical activity and stress may cause an increase in WBC and differential values.
- Pregnancy (final month) and labor may cause increased WBC levels.
- Patients who have had a splenectomy have a persistent, mild elevation of WBC counts.
- ✎ Drugs that may cause *increased* WBC levels include adrenaline, allopurinol, aspirin, chloroform, epinephrine, heparin, quinine, steroids, and triamterene.
- ✎ Drugs that may cause *decreased* WBC levels include antibiotics, anticonvulsants, antihistamines, antimetabolites, antithyroid drugs, arsenicals, barbiturates, chemotherapeutic agents, diuretics, and sulfonamides.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no
- Blood tube commonly used: lavender

## Abnormal findings

## ▲ Increased WBC count (leukocytosis)

Infection  
Leukemic neoplasia  
Trauma  
Stress  
Tissue necrosis  
Inflammation

## ▼ Decreased WBC count (leukopenia)

Drug toxicity  
Bone marrow failure  
Overwhelming infections  
Dietary deficiency  
Autoimmune disease  
Bone marrow infiltration (e.g., myelofibrosis)  
Congenital marrow aplasia

## ▲ ▼ Increased/decreased differential count

See [Table 39](#).

**TABLE 39** Causes of abnormalities in white blood cell (WBC) and differential counts

Type of WBC	Increased	Decreased
Neutrophils	<p><i>Neutrophilia</i></p> <p>Physical or emotional stress</p> <p>Acute suppurative infection</p> <p>Myelocytic leukemia</p> <p>Trauma</p> <p>Cushing syndrome</p> <p>Inflammatory disorders (e.g., rheumatic fever, thyroiditis, rheumatoid arthritis)</p> <p>Metabolic disorders (e.g., ketoacidosis, gout, eclampsia)</p>	<p><i>Neutropenia</i></p> <p>Aplastic anemia</p> <p>Dietary deficiency</p> <p>Overwhelming bacterial infection (especially in elderly patients)</p> <p>Viral infections (e.g., hepatitis, influenza, measles)</p> <p>Radiation therapy</p> <p>Addison disease</p> <p>Drug therapy: myelotoxic drugs (i.e., chemotherapy)</p>

*Continued*

**TABLE 39** Causes of abnormalities in white blood cell (WBC) and differential counts—cont'd

Type of WBC	Increased	Decreased
Lymphocytes	<i>Lymphocytosis</i>	<i>Lymphocytopenia</i>
	Chronic bacterial infection	Leukemia
	Viral infection (e.g., mumps, rubella)	Sepsis
	Lymphocytic leukemia	Immunodeficiency diseases
	Multiple myeloma	Systemic lupus erythematosus
	Infectious mononucleosis	Later stages of HIV infection
	Radiation	Drug therapy: adrenocorticosteroids, antineoplastics
	Infectious hepatitis	Radiation therapy
	Monocytes	<i>Monocytosis</i>
Chronic inflammatory disorders		Aplastic anemia
Viral infections (e.g., infectious mononucleosis)		Hairy-cell leukemia
Tuberculosis		Drug therapy: prednisone
Chronic ulcerative colitis		
Parasites (e.g., malaria)		
Eosinophils	<i>Eosinophilia</i>	<i>Eosinopenia</i>
	Parasitic infections	Increased adrenosteroid production
	Allergic reactions	
	Eczema	
	Leukemia	
	Autoimmune diseases	
Basophils	<i>Basophilia</i>	<i>Basopenia</i>
	Myeloproliferative disease (e.g., myelofibrosis, polycythemia rubra vera)	Acute allergic reactions
	Leukemia	Hyperthyroidism
	Uremia	Stress reaction

## white blood cell scan (WBC scan, Inflammatory scan)

**Type of test** Nuclear scan

### Normal findings

No signs of WBC localization outside the liver or spleen

### Test explanation and related physiology

This test is based on the fact that WBCs are attracted to an area of infection or inflammation. When a patient is suspected of having had infection or inflammation yet the site cannot be localized, the injection of radiolabeled WBCs may identify and localize the area of inflammation or infection. This is especially helpful in patients who have a fever of unknown origin, suspected occult intraabdominal infection, or suspected (yet radiographically inapparent) osteomyelitis. The scan can differentiate infectious from noninfectious processes. For example, it is used to indicate whether an abnormal mass (e.g., a pancreatic pseudocyst) is infected. Areas of noninfectious inflammation (e.g., inflammatory bowel disease) also take up radiolabeled WBCs.

This scan requires drawing blood from the patient, separating out the WBCs, labeling the WBCs with technetium or indium, and reinjecting them back into the patient. Imaging of the whole body 4 to 24 hours later may show an area of increased radioactivity suggestive of accumulation of the radiolabeled WBCs in an area of infection or inflammation.

### Procedure and patient care

#### Before

- PT** Explain the procedure to the patient. See p. xviii for radiation exposure and risks.
- PT** Assure the patient that he or she will not be exposed to large amounts of radioactivity because only tracer doses of the isotope are used.
- PT** Tell the patient that no preparation or sedation is required.

#### During

- Note the following procedural steps:
  1. Approximately 40 to 50 mL of blood is withdrawn from the patient, and the WBCs are extracted from the rest of the blood cells. The WBCs are suspended in saline and tagged with <sup>99m</sup>technetium (<sup>99m</sup>Tc) or <sup>111</sup>indium (<sup>111</sup>In) lipid-soluble product.
  2. The tagged WBCs are reinjected into the patient.

## 980 white blood cell scan

3. At 4, 24, and 48 hours after injection, a gamma ray detector/camera is placed over the body.
4. The patient is placed in supine, lateral, and prone positions so that all surfaces of the body can be visualized.
5. The radionuclide image is recorded on film.

### After

**PT** Inform the patient that because only tracer doses of radioisotopes are used, no precautions need to be taken against radioactive exposure.

### Abnormal findings

- Infection (e.g., abscess or osteomyelitis)
- Inflammation (e.g., inflammatory bowel disease, arthritis)

### notes



## wound culture and sensitivity (C&S)

**Type of test** Microscopic examination

### Normal findings

Negative

### Test explanation and related physiology

Wound cultures are obtained to determine the presence of pathogens in patients with suspected wound infections.

All cultures should be performed before antibiotic therapy is initiated. Otherwise, the antibiotic may interrupt the growth of the organism in the laboratory. More often than not, however, the physician will want to institute antibiotic therapy before the culture results are reported. In these instances a *Gram stain* of the specimen smeared on a slide is most helpful and can be reported in less than 10 minutes. All forms of bacteria are grossly classified as gram positive (blue staining) or gram negative (red staining). Knowledge of the shape of the organism (e.g., spheric, rod shaped) also may be very helpful in the tentative identification of the infecting organism. With knowledge of the Gram stain results, the physician can institute a reasonable antibiotic regimen based on experience of the organism's possible identity.

Most organisms require approximately 24 hours to grow in the laboratory, and a preliminary report can be given at that time. Usually 48 to 72 hours is required for growth and identification of the organism. Cultures may be repeated after appropriate antibiotic therapy to assess for complete resolution of the infection. The *sensitivity* of a bacterium to antibiotics can be determined by adding small antibiotic infused discs to the culture to identify bacterial growth inhibition. The amount of inhibition can be quantitated and provide an accurate list of appropriate antibiotics to which the bacteria will certainly be sensitive.

### Interfering factors

✘ Drugs that may alter test results include antibiotics.

### Procedure and patient care

#### Before

**PT** Explain the procedure to the patient.

#### During

- Aseptically place a sterile cotton swab into the pus of the patient's wound and then place the swab into a sterile, covered

## 982 wound culture and sensitivity

test tube. (Culturing specimens from the skin edge is much less accurate than culturing the suppurative material.)

- If an anaerobic organism is suspected, obtain an anaerobic culture tube from the microbiology laboratory.
- If wound cultures are to be obtained on a patient requiring wound irrigation, obtain the culture before the wound is irrigated.
- If any antibiotic ointment or solution has been previously applied, remove it with sterile water or saline before obtaining the culture.
- Handle all specimens as though they were capable of transmitting disease.
- Indicate on the laboratory slip any medications the patient may be taking that could affect test results.

### After

- Transport the specimen to the laboratory immediately after testing (at least within 30 minutes).
- Notify the physician of any positive results so that appropriate antibiotic therapy can be initiated.

### Abnormal findings

Wound infection

notes

**D-xylose absorption test (Xylose tolerance test)****Type of test** Blood; urine**Normal findings**

Age	60-min plasma (mg/dL)	120-min plasma (mg/dL)	Urine (g/5 hr) [%]
Child	>15-20	>20	>4 [16-32]
Adult	20-57	30-58	>3.5-4 [>14]

**Test explanation and related physiology**

D-Xylose is a monosaccharide that is easily absorbed by the normal intestine. In patients with malabsorption, intestinal D-xylose absorption is diminished; as a result, blood levels and urine excretion are reduced. D-Xylose is the monosaccharide chosen for the test because it is not metabolized by the body. Its serum levels directly reflect intestinal absorption.

This particular monosaccharide does not require pancreatic or biliary exocrine function. Its absorption is directly determined by the small intestine. This test is used to separate patients with diarrhea caused by maldigestion (pancreatic or biliary dysfunction) from those with diarrhea caused by malabsorption (sprue, Whipple disease, Crohn disease).

In this test, the patient is asked to drink a fluid containing a prescribed amount of D-xylose. Blood and urine levels are subsequently evaluated. Excellent gastrointestinal absorption is documented by high blood levels and good urine secretion of D-xylose. Poor intestinal absorption is marked by decreased blood levels and urine excretion.

**Contraindications**

- Patients with abnormal kidney function
- Patients who are dehydrated

**Interfering factors**

- ✚ Drugs that may affect test results include aspirin, atropine, and indomethacin.

## Procedure and patient care

### Before

- PT** Explain the procedure to the patient.
- PT** Instruct the adult patient to fast for 8 hours before testing.
- PT** Tell the pediatric patient or the parents that the patient should fast for at least 4 hours before testing.

### During

- Collect a venous blood sample in a red-top tube before the patient ingests the D-xylose.
  - Collect a first-voided morning urine specimen and send it to the laboratory.
  - Ask the patient to drink the prescribed dose of D-xylose dissolved in 8 oz of water. Record the time of ingestion.
  - Calibrate pediatric doses according to body weight.
  - Repeat venipunctures to obtain blood in exactly 2 hours for an adult and 1 hour for a child.
  - Collect urine for a designated time, usually 5 hours. Refrigerate the urine during the collection period.
  - Observe the patient for nausea, vomiting, and diarrhea, which may occur as side effects of D-xylose.
- PT** Instruct the patient to remain in a restful position. Intense physical activity may alter digestion and affect results.

### After

- Apply pressure to the venipuncture site.
- PT** Inform the patient that normal activity may be resumed after completion of the study.

## Abnormal findings

### ▼ Decreased levels

- Sprue
- Lymphatic obstruction
- Enteropathy (e.g., radiation)
- Crohn disease
- Whipple disease
- Small intestine bacterial overgrowth
- Hookworm
- Viral gastroenteritis
- Giardia lamblia infestation
- Short-bowel syndrome

## Zika virus (ZIKV)

**Type of test** Blood, urine, cerebrospinal fluid (CSF)

### Normal findings

No Zika detected

### Test explanation and related physiology

The Zika virus is transmitted to humans by an infected *Aedes* species of mosquito. Sexual transmission among humans has also been described. This testing is used to identify Zika infections in people who have symptoms of the infection and live in or who have traveled to an area where Zika is known to exist. A blood or urine test can be used to confirm Zika virus infection.

Zika virus belongs to the genus *Flavivirus* single-stranded RNA virus. In most cases, Zika virus infection causes a mild, self-limited illness. Symptoms include fever, arthralgia, retroocular headache, and conjunctivitis. Typically, the predominant symptom is a rash. Symptoms last from 2 to 7 days. Although Zika virus infection is generally well tolerated, infection in pregnant women can lead to microcephaly and other neurologic defects in fetuses because of transplacental transmission. Infection may precipitate a decision for an early abortion.

Diagnosis of Zika virus infection is typically based on serologic tests that identify immunoglobulin M (IgM) antibodies to Zika. IgM occurs early in the clinical course and can be found in the serum and CSF. Because viral levels may be higher in urine and last longer than in serum, urine testing for Zika has been used.

Molecular testing (reverse-transcriptase polymerase chain reaction) of serum and urine can detect Zika RNA early in the course of disease. Zika viral IgM testing (enzyme-linked immunosorbent assay) can provide presumptive confirmation of Zika disease. Because symptoms of Zika are so similar to other viral diseases, multiviral panels have been developed to test for other diseases, such as Dengue and Chikungunya.

### Interfering factors

- Serologic testing for Zika can cross react with other flavivirus infections, such as yellow fever, West Nile virus, and dengue.

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: no

## 986 Zika virus

- Blood tube commonly used: lavender and red
- See inside front cover for Routine Urine Testing.
- Transfer urine to a clean vial with a screw cap and O-ring to prevent leakage.

### **Abnormal findings**

Zika infection

notes

## zinc protoporphyrin (ZPP)

**Type of test** Blood

### Normal findings

0-69  $\mu\text{mol ZPP/mol heme}$

### Test explanation and related physiology

ZPP is used in screening for iron deficiency anemia or lead poisoning. It is also used in monitoring treatment of chronic lead poisoning. ZPP is found in red blood cells when heme production is inhibited by lead toxicity. Lead prevents iron, but not zinc, from attaching to the protoporphyrin. If there is iron deficiency, instead of incorporating a ferrous ion to form heme, protoporphyrin (the immediate precursor of heme) incorporates a zinc ion, forming ZPP. In addition to lead poisoning and iron deficiency, zinc protoporphyrin levels can be elevated as the result of a number of other conditions (e.g., sickle cell anemia).

### Procedure and patient care

- See inside front cover for Routine Blood Testing.
- Fasting: yes (12 hours)
- Blood tube commonly used: verify with laboratory.

### Abnormal findings

#### ▲ Increased levels

Lead poisoning  
Iron deficiency  
Anemia of chronic illness  
Sickle cell anemia  
Sideroblastic anemia  
Vanadium exposure

notes

## appendix A: list of tests by body system

Tests in this list are grouped by the following: cancer studies; cardiovascular, endocrine, gastrointestinal, hematologic, hepatobiliary, and immunologic systems; miscellaneous studies; and nervous, pulmonary, renal/urologic, reproductive, and skeletal systems.

### CANCER STUDIES

- Acid phosphatase, 7
- Bence Jones protein, 134
- Beta<sub>2</sub>-microglobulin, 628
- Bladder cancer markers, 145
- Bone scan, 165
- Breast cancer genetic testing, 445
- Breast cancer genomics, 454
- Breast cancer tumor analysis, 176
- Breast ductal lavage, 179
- CA 15-3 and CA 27.29 tumor markers, 194
- CA 19-9 tumor marker, 194
- CA125 tumor marker, 194
- Carcinoembryonic antigen, 194
- Cathepsin D, 176
- Cell culture drug resistance testing, 216
- Cervical biopsy, 224
- Colon cancer genetic testing, 445
- Colon cancer tumor analysis, 256
- Des gamma carboxy prothrombin, 195
- DNA ploidy status, 176
- Ductoscopy, 336
- Early prostate cancer antigen, 743
- Estrogen receptor assay, 389
- Gallium scan, 433
- HER 2 protein, 176
- Ki67 protein, 176
- Liquid biopsy, 569
- Mammography, 606
- Melanoma genetic testing, 445
- Microglobulin, 628
- Neuron specific enolase, 638
- Octreotide scan, 652
- Ovarian cancer genetic testing, 445
- p53 protein, 176
- Papanicolaou smear, 667
- Progesterone receptor assay, 735
- ProstaScint scan, 739
- Prostate-specific antigen, 743
- Salivary gland nuclear imaging, 795
- Sentinel lymph node biopsy, 805
- Serotonin, 807
- S-phase fraction, 176
- Sputum cytology, 843
- Thyroid cancer genetic testing, 445
- Thyroid cancer genomic testing, 453
- Thyroid fine needle aspiration biopsy, 881
- Tumor markers, 194



**CARDIOVASCULAR SYSTEM**

- Adenosine stress test, 207  
 Aldosterone, 25  
 Antimyocardial antibody, 83  
 Antistreptolysin O titer, 850  
 Apolipoproteins, 101  
 Arteriography, 114  
 Aspartate aminotransferase, 125  
 Atrial natriuretic peptide, 636  
 Brain natriuretic peptide, 636  
 Cardiac catheterization, 201  
 Cardiac MRI, 600  
 Cardiac nuclear scanning, 207  
 Cardiac stress testing, 210  
 Cardiac genetic testing, 450  
 Carotid artery duplex scanning, 214  
 Catecholamines, 954  
 Ceramides, 222  
 Chest x-ray, 227  
 CHF peptides, 636  
 Cholesterol, 235  
 Computed tomography of the chest, 281  
 Computed tomography of the heart, 284  
 Creatine kinase, 297  
 Creatine phosphokinase, 297  
 Cryoglobulin, 306  
 C-type natriuretic peptide, 636  
 Digital subtraction angiography, 114  
 Dipyridamole stress test, 211  
 Dobutamine stress test, 211  
 Echocardiography, 339  
 Electrocardiography, 342  
 Electrophysiologic study, 359  
 Fibrinogen, 422  
 Galectin-3, 430  
 Holter monitoring, 510  
 Homocysteine, 513  
 Intravascular ultrasound, 543  
 Ischemia-modified albumin, 550  
 Isonitrile scan, 207  
 Lactic dehydrogenase, 553  
 Lipoprotein-associated phospholipase A<sub>2</sub>, 564  
 Lipoproteins, 565  
 Microvolt T-wave alternans, 344  
 MUGA scan, 207  
 Myoglobin, 635  
 Natriuretic peptide, 636  
 N-terminal fragment of pro-brain natriuretic peptide, 636  
 Pericardiocentesis, 689  
 Plethysmography, arterial, 713  
 Positron emission tomography, 719  
 Renin assay, 783  
     Plasma, 783  
     Renal vein, 783  
 Tilt-table test, 359  
 Transesophageal echocardiography, 903  
 Transthoracic echocardiography, 339  
 Triglycerides, 908  
 Troponins, 912  
 Vanillylmandelic acid and catecholamines, 954  
 Vascular ultrasound studies, 960  
 Venography of lower extremities, 963

## ENDOCRINE SYSTEM

- Adrenocorticotrophic hormone
  - Stimulation test, 13
  - With cosyntropin, 13
  - With metyrapone, 16
- Aldosterone, 25
- Androstenediones, 18
- Antidiuretic hormone, 71
- Antithyroglobulin antibody, 99
- Antithyroid peroxidase
  - antibody, 100
- Blood glucose, 462
- Calcitonin, 187
- Calcium, 189
- Catecholamines, 954
- Chromosome karyotype, 244
- Computed tomography of
  - adrenals, 272
- Cortisol, blood and urine, 290
- C-peptide, 293
- Dehydroepiandrosterone, 18
- Dehydroepiandrosterone
  - sulfate, 18
- Dexamethasone suppression
  - test, 324
- Diabetes mellitus
  - autoantibody panel, 327
- Erythropoietin, 375
- Estriol excretion, 386
- Estrogen fractions, 386
- Follicle-stimulating hormone
  - assay, 593
- Gastrin, 439
- Glucagon, 460
- Glucose
  - Blood, 462
  - Mean plasma, 472
  - Postprandial, 465
  - Tolerance test, 467
  - Urine, 467
- Glutamic acid decarboxylase
  - antibody, 327
- Glycated proteins, 471
- Glycosylated hemoglobin, 471
- Growth hormone, 474
  - Stimulation test, 477
  - Suppression test, 475
- 17-Hydroxycorticosteroids, 526
- 21-Hydroxylase antibodies, 530
  - Insulin assay, 539
  - Insulin autoantibody, 327
  - Insulin-like growth factor
    - binding proteins, 541
  - Islet cell antibody, 327
- Ketones, 930
- 17-Ketosteroids, 862
- Long-acting thyroid
  - stimulator, 889
- Luteinizing hormone assay, 593
- Metanephrine, plasma free, 621
- Metyrapone test, 16
- Osmolality, blood, 654
- O'Sullivan test, 465
- Parathyroid, 675
  - Hormone, 675
  - Scan, 677
- Pheochromocytoma
  - suppression and provocative
    - testing, 692
- Phosphate, 694
- Phosphorus, 694
- Postprandial glucose, 465
- Prolactin levels, 737
- Renin assay, 783
  - Plasma, 783
  - Renal vein, 784
- Reverse triiodothyronine
  - (T<sub>3</sub>), 910
- Somatomedin C, 541
- Testosterone, 861
- Thyroglobulin, 879

- Thyroid-binding inhibitory immunoglobulins, 889  
 Thyroid cancer genetic testing, 445  
 Thyroid cancer genomic testing, 453  
 Thyroid fine needle aspiration biopsy, 881  
 Thyroid scanning, 883  
 Thyroid-stimulating hormone, 886  
 Thyroid-stimulating immunoglobulins, 889  
 Thyroid ultrasound, 891  
 Thyrotropin receptor antibody, 889  
 Thyrotropin-releasing hormone stimulation test, 892  
 Thyroxine ( $T_4$ ), free and total, 894  
 Thyroxine-binding globulin, 897  
 TRH stimulation test, 892  
 Triiodothyronine ( $T_3$ ), 910  
 TSH stimulation test, 888  
 Vanillylmandelic acid (VMA), 954

### GASTROINTESTINAL SYSTEM

- Abdominal ultrasound, 1  
 Antiparietal cell antibody, 90  
 Barium  
     Enema, 128  
     Swallow, 131  
 Capsule endoscopy, 382  
 Carcinoembryonic antigen, 194  
*Clostridium difficile* testing, 246  
 Colon cancer tumor analysis, 256  
 Colonoscopy, 258  
 Computed tomography  
     Arteriography, 284  
     Colonography, 273  
     Of the abdomen, 272  
 Cytolethal distending toxin B, 316  
 DNA stool sample, 844  
 Endomysial antibodies, 458  
 Esophageal electrical impedance, 378  
 Esophageal function studies, 377  
 Esophagogastroduodenoscopy, 382  
 Fecal fat, 402  
 Fecal calprotectin, 400  
 Fecal immunochemical test, 844  
 Gastric emptying scan, 437  
 Gastrin, 439  
 Gastroesophageal reflux scan, 441  
 Gastrointestinal bleeding scan, 443  
 Gliadin antibodies, 458  
*Helicobacter pylori* antibodies test, 482  
 5-Hydroxyindoleacetic acid, 528  
 Lactoferrin, 555  
 Lactose  
     Breath test, 557  
     Tolerance test, 557  
 Meckel diverticulum nuclear scan, 616  
 Obstruction series, 650  
 Paracentesis, 671  
 Pepsinogen, 688  
 Prealbumin, 728  
 Protein, blood, 746  
 Septin 9 DNA methylation assay, 625

## 992 appendix A: list of tests by body system

- Sialography, 818
- Sigmoidoscopy, 258
- Small bowel follow-through, 830
- Small intestinal bacterial overgrowth, 833
- Stool culture, 847
- Stool for leukocytes, 849
- Stool for occult blood, 844
- Swallowing examination, 855
- Tissue transglutaminase antibodies, 458
- Transrectal ultrasonography, 741
- Upper gastrointestinal x-ray study, 922
- Urea breath test, 925
- Urea nitrogen blood test, 155
- Videofluoroscopy, 855
- Virtual colonoscopy, 258
- D-Xylose absorption test, 985

### HEMATOLOGIC SYSTEM

- Anti-Factor Xa, 77
- Activated clotting time, 9
- Activated protein C resistance, 681
- Antithrombin III, 97
- Antithrombin activity and antigen assay, 97
- Basophils, 974
- Blood smear, 149
- Blood typing, 152
- Bone marrow biopsy, 160
- Coagulating factors concentration, 248
- Complete blood count and differential count, 269
- Coombs test,
  - Direct, 287
  - Indirect, 289
- D-Dimer test, 319
- Delta-aminolevulinic acid, 321
- 2,3-Diphosphoglycerate, 328
- Disseminated intravascular coagulation screening, 330
- Eosinophils, 974
- Erythrocyte fragility, 371
- Erythropoietin, 375
- Factor V Leiden, 397
- Ferritin, 404
- Fibrin monomers, 877
- Fibrinogen, 422
- Fibrinopeptide A, 877
- Folic acid, 426
- Haptoglobin, 479
- Hematocrit, 485
- Hemochromatosis genetic testing, 445
- Hemoglobin, 488
  - Electrophoresis, 490
- Intrinsic factor antibody, 545
- Iron level and total iron-binding capacity, 546
- Leukoagglutinin test, 639
- Lymphocytes, 974
- Methemoglobin, 623
- Monocytes, 974
- Neutrophil antibody screen, 639
- Neutrophils, 974
- Partial thromboplastin time, activated, 681
- PI-linked antigen, 697
- Plasminogen, 699
- Plasminogen activator inhibitor I, 700
- Platelet, 702
  - Aggregation test, 702
  - Antibody detection, 704
  - Closure time, 709
  - Count, 706
  - Function assay, 709
  - Volume, mean, 712

- Porphyrins and porphobilinogens, 717
- Protein C, 751
- Protein S, 751
- Prothrombin, 753
  - Fragment, 877
  - Time, 753
- Red blood cell, 770
  - Casts, urine, 936
  - Count, 770
  - Indices, 772
- Reticulocyte count, 788
- Ristocetin test, 702
- Sickle cell screen, 820
- Thromboelastography, 874
- Thrombosis indicators, 877
- Total blood volume, 900
- Transferrin receptor assay, 906
- Urinary methylmalonic acid, 969
- Uroporphyrinogen-1-synthase, 953
- Vitamin B<sub>12</sub>, 969
- Warfarin metabolism genetic test panel, 755
- White blood cell count and differential count, 974
- Zinc protoporphyrin, 987

### HEPATOBIILIARY SYSTEM

- Abdominal ultrasound, 1
- Alanine aminotransferase, 21
- Aldolase, 23
- Alkaline phosphatase, 29
- Alpha-fetoprotein, 41
- Ammonia level, 47
- Amylase
  - Blood, 56
  - Urine, 56
- Anti-liver/kidney microsomal type 1 antibodies, 81
- Aspartate aminotransferase, 125
- Bilirubin, 137
- CA 19-9 tumor marker, 194
- Cholesterol, 235
- Computed tomography of the abdomen, 272
- Endoscopic retrograde cholangiopancreatography, 364
- Epstein-Barr virus titer, 368
- Fecal fat, 402
- Gallbladder nuclear scanning, 431
- Gallium scan, 433
- Gamma-glutamyl transpeptidase, 435
- Heinz body preparation, 481
- Hepatitis B DNA testing, 493
- Hepatitis C RNA testing, 493
- Hepatitis D viral antigen testing, 493
- Hepatitis virus studies, 493
- Lactic dehydrogenase, 553
- Leucine aminopeptidase, 561
- Lipase, 562
- Liver biopsy, 571
- Liver and pancreaticobiliary system ultrasonography, 1
- Liver/spleen scanning, 574
- 5'-Nucleotidase, 649
- Obstruction series, 650
- Pancreatic elastase, 660
- Pancreatic enzymes, 662
- Pancreatobiliary FISH testing, 665
- Percutaneous transhepatic cholangiography, 364
- Protein, 746
- Secretin-pancreozymin, 662
- Sweat electrolytes test, 857

## IMMUNOLOGIC SYSTEM

- Agglutinins,
  - Cold, 253
  - Febrile, 399
- Aldolase, 23
- Allergy blood testing, 31
- Allergy skin testing, 33
- Anticardiolipin antibodies, 65
- Anticentromere antibody test, 67
- Antichromatin antibody, 68
- Anti-cyclic citrullinated peptide antibody, 70
- Antideoxyribonuclease-B titer, 850
- Anti-DNA antibody test, 74
- Anti-extractable nuclear antigens, 75
- Anti-glomerular basement membrane antibodies, 79
- Anti-glycan antibodies, 80
- Antihistone antibody, 68
- Anti-Jo-1 antibodies, 75
- Anti-liver/kidney type 1 antibodies, 81
- Antimitochondrial antibody, 82
- Antimyocardial antibody, 83
- Antineutrophil cytoplasmic antibody, 84
- Antinuclear antibody, 86
- Antinucleosome antibody, 68
- Antiparietal cell antibody, 90
- Antiscleroderma antibody, 91
- Antismooth muscle antibody, 92
- Anti-SS-A, anti-SS-B, and anti-SS-C antibodies, 95
- 11 Beta-prostaglandin F(2) alpha, 136
- CD4/CD8 for HIV, 219
- Cell surface immunophenotyping, 219
- Cold agglutinins, 253
- C-reactive protein, 295
- Complement assay, 266
- Cryoglobulin, 306
- Cutaneous immunofluorescence biopsy, 824
- Cytokines, 314
- Diabetes mellitus autoantibody panel, 327
- Epstein-Barr virus titer, 368
- Febrile agglutinins, 399
- Fungal testing, 428
- Glutamic decarboxylase antibody, 327
- Heparin-induced thrombocytopenia antibodies, 704
- HIV
  - Drug resistance testing, 502
  - Oral testing, 508
  - Serology, 506
  - Urine testing, 508
  - Viral load, 503
- Human lymphocyte antigen B27, 518
- Human T-cell lymphotropic virus I/II antibody, 525
- Immunofixation electrophoresis, 536
- Insulin autoantibody, 327
- Islet cell antibody, 327
- Lyme disease test, 595
- Measles rubeola antibody, 615
- Microglobulin, 628
- Mononucleosis spot test, 630
- Parvovirus B19 antibody, 684
- Protein electrophoresis, 746
- Rabies-neutralizing antibody test, 769

Rheumatoid factor, 790	<i>Streptococcus</i> serologic testing, 850
Ribosome P antibodies, 792	<i>Varicella</i> virus, 958
Rubella antibody, 793	Virus testing, 967
Rubeola antibody, 615	Zika virus, 985
Skin biopsy, 824	

### MISCELLANEOUS STUDIES

Age-related macular degeneration risk analysis, 20	Nicotine and metabolites, 646
Aluminum, 43	Paternity genetic testing, 445
Anion gap, 63	PET/CT image fusion, 174
Bioterrorism infectious agent testing, 141	Positron emission tomography (PET), 719
Blood culture and sensitivity, 147	Potassium, 724
Carbon dioxide content, 197	Blood, 724
Carboxyhemoglobin, 199	Urine, 727
Chloride, blood, 233	Protein, 746
Cholinesterase, 238	Rubeola antibody, 615
Cotinine, 646	Sexual assault testing, 809
C-reactive protein test, 295	Sialography, 818
Cytochrome P450 genotyping, 334	Sleep studies, 826
Dental x-rays, 322	Sodium, 835
Drug monitoring, 332	Blood, 835
Drug sensitivity genotype, 332	Fractional excretion, 837
Erythrocyte sedimentation rate, 373	Urine, 837
Ethanol, 391	Substance abuse testing, 852
Fluorescein angiography, 424	Tay-Sachs disease genetic testing, 445
Forensic genetic testing, 445	Therapeutic drug monitoring, 332
Fractional excretion of sodium, 837	Throat and nose cultures, 872
Fungal antibody tests, 428	Toxicology screening, 852
Genetic testing, 445	Urine culture and sensitivity, 947
Genomic testing, 453	Varicella virus, 958
Glucose	Viral cultures, 965
Blood, 462	Virus testing, 967
Urine, 467	Virtual autopsy, 274
Lead, 620	West Nile virus testing, 967
Magnesium, 597	White blood cell scan, 979
Magnetic resonance testing (MRI), 599	Wound culture and sensitivity, 981
	Zika virus, 985

## NERVOUS SYSTEM

- Acetylcholine receptor antibody, 5
- Amyloid beta protein precursor, soluble, 58
- Apolipoproteins, 101
- Apolipoprotein E-4, 102
- Brain scan, 174
- Caloric study, 192
- Cerebrospinal fluid examination, 576
- Computed tomography of the brain, 278
- C-reactive protein test, 295
- Digital subtraction angiography, 114
- Electrocorticography, 347
- Electroencephalography, 347
- Electromyography, 350
- Electroneurography, 355
- Electronystagmography, 357
- Evoked potential studies, 393
- Helical CT scan, brain, 278
- Hexosaminidase, 500
- Lumbar puncture, 576
- Magnetoencephalography, 347
- Magnetic resonance angiography, 600
- Magnetic resonance imaging, 599
- Magnetic resonance spectroscopy, 599
- Magnetic resonance venography, 601
- Myelography, 632
- Parkinson disease testing, 679
- Positron emission tomography (PET), 719
- Skull x-ray, 825
- Spinal x-ray, 839
- Tau protein, 58
- Tay-Sachs disease genetic testing, 445

## PULMONARY SYSTEM

- Acid-fast bacillus smear, 917
- Alpha<sub>1</sub>-antitrypsin, 37
- Alpha<sub>1</sub>-antitrypsin phenotyping, 37
- Angiotensin-converting enzyme, 61
- Arterial blood gases, 106
- Blood gases, 106
- Body plethysmography, 762
- Bronchoscopy, 183
- Carbon dioxide content, 197
- Carboxyhemoglobin, 199
- Chest x-ray, 227
- Computed tomography of the chest, 281
- Cystic fibrosis genetic testing, 445
- Gas dilution test, 762
- Laryngoscopy, 184
- Legionnaires disease antibody test, 560
- Lung biopsy, 584
- Lung cancer molecular testing, 588
- Lung scan, 590
- Mediastinoscopy, 618
- Mycoplasma pneumoniae* antibodies, 631
- Nicotine and metabolites, 646
- Nose culture, 872
- Oximetry, 658
- Pleural biopsy, 715
- Pulmonary angiography, 757
- Pulmonary function tests, 759



QuantiFERON-TB Gold test, 919	Thoracentesis and pleural fluid analysis, 864
SARS testing, 797	Thoracoscopy, 870
Sleep studies, 826	Throat and nose cultures, 872
Sputum culture and sensitivity, 841	Tuberculin test, 915
Sputum cytology, 843	Tuberculosis culture, 917
<i>Streptococcus</i> screen, 850	Tuberculosis testing, 919
	Urine pH, 927

### RENAL/UROLOGIC SYSTEM

Acid phosphatase, 7	Fluorescence in situ hybridization (FISH) bladder molecular genetic testing, 145
Aldosterone, 25	Granular casts, 936
Amino acid profiles, 45	Hyaline casts, 936
Angiotensin, 60	Intravenous pyelography, 764
Anion gap, 63	Ketones, 935
Antidiuretic hormone, 71	Kidney sonogram, 1
Antistreptolysin O titer, 850	Kidney, ureter, bladder x-ray study, 650
Beta <sub>2</sub> -microglobulin, 628	Lactic acid, 551
Bilirubin, 137	Leukocyte esterase, 930
Bladder cancer markers, 145	Microalbumin, 626
Bladder tumor antigen, 145	Neutrophil gelatinase-associated lipocalin, 640
Captopril renal scan, 779	Nitrites, 935
Carbon dioxide content, 197	Nuclear matrix protein, 145
Catecholamines, 954	Osmolality
Chloride, blood, 233	Blood, 654
Computed tomography of kidney, 274	Urine, 656
Creatinine, blood, 301	Pelvic floor sphincter electromyography, 353
Creatinine clearance, 303	Phospholipase A2 receptor antibody, 696
Crystals, 936, 939	Potassium, 724
Cystatin C, 302	Blood, 724
Cystography, 307	Urine, 727
Cystometry, 949	ProstaScint scan, 739
Cystoscopy, 309	Prostate/rectal sonogram, 741
Endourethral urologic ultrasound, 1	Prostate-specific antigen, 743
Epithelial casts, 937	
Erythropoietin, 375	
Estimated glomerular filtration rate, 303	
Fatty casts, 936	

## 998 appendix A: list of tests by body system

- Protein, urine, 746
- PSA velocity, 743
- Pyelography, 764
- Renal angiography, 114
- Renal biopsy, 776
- Renal scanning, 779
- Renin assay, 783
  - Plasma, 783
  - Renal vein, 783
- Retrograde pyelography, 764
- Scrotal ultrasound, 799
- Sodium, 835
  - Blood, 835
  - Urine, 837
- Testicular ultrasound, 799
- Tubular casts, 937
- Urea nitrogen blood test, 155
- Urethral pressure profile, 949
- Uric acid, 927
- Urinalysis, 930
- Urinary stone analysis, 945
- Urine, 927
  - Appearance and color, 930
  - Culture and sensitivity, 947
  - Flow studies, 949
  - Odor, 930
  - pH, 930
  - Specific gravity, 930
- Urodynamic studies, 949
- Vanillylmandelic acid, 954
- Waxy casts, 936
- White blood cells and casts, urine, 930

### REPRODUCTIVE SYSTEM

- Alpha-fetoprotein, 41
- Amniocentesis, 49
- Amniotic fluid index, 407
- Antispermatozoal antibody, 93
- Apt test, 105
- Biophysical profile, fetal, 406
- Breast cancer
  - Genetic testing, 445
  - Tumor analysis, 178
- Breast sonogram, 181
- Breast ultrasonography, 181
- CA 15-3 and CA 27.29 tumor markers, 194
- CA125 tumor marker, 194
- Cell-free maternal DNA testing, 217
- Cervical biopsy, 224
- Chlamydia*, 230
- Chorionic villus sampling, 241
- Colposcopy, 262
- Contraction stress test, 409
- Cytomegalovirus, 317
- Electromyography of the pelvic floor sphincter, 353
- Endometrial biopsy, 362
- Estriol excretion, 386
- Estrogen fractions, 386
- Fetal
  - Biophysical profile, 406
  - Contraction stress test, 409
  - Fibronectin, 412
  - Hemoglobin, 413
  - Nonstress test, 415
  - Nuchal translucency, 41
  - Oxygen saturation monitoring, 417
  - Scalp blood pH, 417
- Fetoscopy, 419
- Follicle-stimulating hormone, 593
- Gonorrhea culture, 813
- Herpes genitalis, 498
- Human chorionic gonadotropin, 515
- Human papillomavirus, 520
- Human placental lactogen, 523
- Hyaluronan binding assay, 803

- Hysterosalpingography, 531
- Hysteroscopy, 533
- Inhibin A, 613
- IUD localization, 534
- Lamellar body count, 49
- Laparoscopy, 535
- Luteinizing hormone assay, 593
- Mammography, 606
- Maternal screen testing, 612
- Newborn metabolic screening, 642
- Noninvasive prenatal testing, 217
- Nonstress test fetal, 415
- Obstetric ultrasonography, 685
- Oxytocin challenge test, 409
- Papanicolaou smear, 667
- Pelvic ultrasonography, 685
- Phenylketonuria test, 642
- Placental alpha microglobulin, 52
- Placental growth factor, 698
- Pregnancy-associated plasma protein-A, 730
- Pregnancy tests, 515
- Pregnanediol, 731
- Progesterone assay, 733
- Prolactin levels, 737
- Rubella antibody test, 793
- Saline infusion sonography, 686
- Semen analysis, 801
- Sexual assault testing, 809
- Sexually transmitted disease testing, 813
- Sims-Huhner test, 822
- Sperm chromatin structure assay, 802
- Sperm DNA fragmentation assay, 802
- Sperm DNA integrity, 802
- Sperm functions test, 802
- Sperm penetration assay, 802
- Syphilis detection test, 859
- TORCH test, 899
- Toxoplasmosis antibody titer, 902

### SKELETAL SYSTEM

- Aldolase, 23
- Alkaline phosphatase, 29
- Alpha defensin test, 39
- Arthrocentesis with synovial fluid analysis, 118
- Arthroscopy, 122
- Bone
  - Densitometry, 157
  - Scan, 165
  - Turnover markers, 168
  - X-ray, 172
- Electromyography, 350
- Human lymphocyte antigen B27, 518
- Magnetic resonance imaging, 599
- Myelography, 632
- Myoglobin, 635
- N-telopeptide, 168
- Osteocalcin, 168
- Pyridinium crosslinks, 168
- Rheumatoid factor, 790
- Spinal x-rays, 839
- Uric acid, 927
- Vitamin D, 971

## appendix B: list of tests by type

Tests in this list are grouped by the following types: blood, electrodiagnostic, endoscopy, fluid analysis, manometric, microscopic examinations, nuclear scans, other studies, sputum, stool, ultrasound, urine, and x-ray.

### BLOOD TESTS

- |   |   |
|---|---|
| Acetylcholine receptor antibody, 5              | Anion gap, 63                                     |
| Acid phosphatase, 7                             | Anticardiolipin antibodies, 65                    |
| Activated clotting time, 9                      | Anticentromere antibody test, 67                  |
| Activated protein C resistance, 397             | Antichromatin antibody test, 68                   |
| Adrenocorticotrophic hormone, 11                | Anti-cyclic citrullinated peptide antibody, 70    |
| Stimulation test with cosyntropin, 13           | Antideoxyribonuclease-B titer, 74                 |
| Stimulation test with metyrapone, 16            | Antidiuretic hormone, 71                          |
| Age-related macular degeneration, 20            | Anti-DNA antibody test, 74                        |
| Agglutinins                                     | Anti-extractable nuclear antigens, 75             |
| Cold, 253                                       | Antifactor Xa, 77                                 |
| Febrile, 399                                    | Anti-glomerular basement membrane antibodies, 79  |
| Alanine aminotransferase, 21                    | Antiglycan antibodies, 80                         |
| Albumin, 746                                    | Antihistone antibodies, 68                        |
| Aldolase, 23                                    | Anti-Jo-1 antibodies, 75                          |
| Aldosterone, 25                                 | Antiliver/kidney microsomal type 1 antibodies, 81 |
| Alkaline phosphatase, 29                        | Antimitochondrial antibody, 82                    |
| Allergy blood testing, 31                       | Antimyocardial antibody, 83                       |
| Alpha <sub>1</sub> -antitrypsin, 37             | Antineutrophil cytoplasmic antibody, 84           |
| Alpha <sub>1</sub> -antitrypsin phenotyping, 37 | Antinuclear antibody, 86                          |
| Alpha-fetoprotein, 41                           | Antinucleosome antibodies, 68                     |
| Aluminum, 43                                    | Antiparietal cell antibody, 90                    |
| Amino acid profiles, 45                         | Antiscleroderma antibody, 91                      |
| Ammonia level, 47                               | Anti-smooth muscle antibody, 92                   |
| Amylase, 56                                     | Antispermatozoal antibody, 93                     |
| Androstenediones, 18                            |   |
| Angiotensin, 60                                 |   |
| Angiotensin-converting enzyme, 61               |   |

- Anti-SS-A, anti-SS-B, and anti-SS-C antibodies, 95
- Antistreptolysin O titer, 850
- Antithrombin III, 97
- Antithrombin activity and antigen assay, 97
- Antithyroglobulin antibody, 99
- Antithyroid peroxidase antibody, 100
- Apolipoproteins, 101
- Arterial blood gases, 106
- Aspartate aminotransferase, 125
- Atrial natriuretic peptide, 636
- Basic metabolic panel (BMP), 133
- Basophils, 974
- Beta<sub>2</sub>-microglobulin, 628
- Bilirubin, 137
- Bioterrorism infectious agents, 141
- Blood
  - Culture and sensitivity, 147
  - Gases, 106
  - Smear, 149
  - Typing, 152
- Blood urea nitrogen (BUN), 155
- Bone turnover markers, 168
- Brain natriuretic peptide, 636
- Breast cancer genetic testing, 445
- CA 15-3 and CA 27.29 tumor markers, 194
- CA 19-9 tumor marker, 194
- CA125 tumor marker, 194
- Calcitonin, 187
- Calcium, 189
- Carbon dioxide content, 197
- Carboxyhemoglobin, 199
- Carcinoembryonic antigen, 194
- CD4/CD8 for HIV, 219
- Cell-free maternal DNA testing, 217
- Cell surface immunophenotyping, 219
- Ceramides, 222
- Chlamydia*, 230
- Chloride, 233
- Cholesterol, 235
- Cholinesterase, 238
- Chromosome karyotype, 244
- Clonidine suppression test, 692
- Coagulating factors
  - concentration, 248
- Coagulation profile, 255
- Cold agglutinins, 253
- Complement assay, 266
- Complete blood count, 269
- Comprehensive metabolic profile, 270
- Coombs test
  - Direct, 287
  - Indirect, 289
- Cortisol, 290
- Cotinine, 646
- C-peptide, 293
- C-reactive protein test, 295
- Creatine kinase, 297
- Creatine phosphokinase, 297
- Creatinine, 301
- Creatinine clearance, 303
- Cryoglobulin, 306
- C-type natriuretic peptide, 636
- Cytokines, 314
- Cytolethal distending toxin B, 316
- Cytomegalovirus, 317
- Dehydroepiandrosterone, 18
- Dehydroepiandrosterone sulfate, 18
- Des gamma carboxy prothrombin, 195
- Dexamethasone suppression test, 324
- D-Dimer test, 319

## 1002 appendix B: list of tests by type

- Diabetes mellitus
  - autoantibody panel, 327
- Differential count, 974
- 2,3-Diphosphoglycerate, 328
- Disseminated intravascular coagulation screening, 330
- Drug monitoring, 332
- Drug sensitivity genotype, 332
- Early prostate cancer antigen, 743
- Endomysial antibodies, 458
- Eosinophils, 974
- Epstein-Barr virus titer, 368
- Erythrocyte
  - Fragility, 371
  - Sedimentation rate, 373
- Erythropoietin, 375
- Estriol excretion, 386
- Estrogen fractions, 386
- Ethanol, 391
- Factor V Leiden, 397
- Febrile agglutinins, 399
- Ferritin, 404
- Fetal hemoglobin, 413
- Fetal scalp blood pH, 417
- Fibrin monomers, 877
- Fibrinogen, 422
- Fibrinopeptide A, 877
- Folic acid, 426
- Follicle-stimulating hormone assay, 593
- Free thyroxine, 894
- Fungal antibody tests, 428
- Galectin-3, 430
- Gamma-glutamyl transpeptidase, 435
- Gastrin, 439
- Genetic testing, 445
- Genomic prostate testing, 454
- Gliadin antibodies, 458
- Glucagon, 460
- Glucagon stimulation test, 692
- Glucose
  - Blood, 462
  - Mean plasma, 472
  - Postprandial, 465
  - Tolerance test, 467
- Glutamic acid decarboxylase antibody, 327
- Glycated proteins, 471
- Glycosylated hemoglobin, 471
- Growth hormone, 474
  - Stimulation test, 477
  - Suppression test, 476
- Haptoglobin, 479
- Heinz body preparation, 481
- Helicobacter pylori* antibody, 482
- Hematocrit, 485
- Hemochromatosis genetic testing, 445
- Hemoglobin, 488
  - Electrophoresis, 490
- Heparin-induced thrombocytopenia antibodies, 704
- Hepatitis B DNA testing, 493
- Hepatitis C RNA testing, 493
- Hepatitis D viral antigen, 493
- Hepatitis virus studies, 493
- Hexosaminidase, 500
- HIV
  - Drug resistance testing, 502
  - Serology, 506
  - Viral load, 503
- HLA-B27 antigen, 518
- Homocysteine, 513
- Human chorionic gonadotropin, 515
- Human lymphocyte antigen B27, 518
- Human placental lactogen, 523
- Human T-cell lymphotropic virus I/II antibody, 525

- 21-Hydroxylase antibodies, 530  
 Inhibin A, 613  
 Insulin assay, 539  
 Insulin autoantibody, 327  
 Insulin-like growth factor binding proteins, 541  
 Intrinsic factor antibody, 545  
 Iron level and total iron-binding capacity, 546  
 Ischemia-modified albumin, 550  
 Islet cell antibody, 327  
 Lactic acid, 551  
 Lactic dehydrogenase, 553  
 Lactose tolerance test, 557  
 Lead, 620  
 Legionnaires disease antibody test, 560  
 Leucine aminopeptidase, 561  
 Leukoagglutinin test, 639  
 Lipase, 562  
 Lipoprotein-associated phospholipase A<sub>2</sub>, 564  
 Lipoproteins, 565  
 Liquid biopsy, 569  
 Long-acting thyroid stimulator, 889  
 Luteinizing hormone assay, 593  
 Lyme disease test, 595  
 Lymphocytes, 974  
 Magnesium, 597  
 Maternal screen testing, 612  
 Measles rubeola antibody, 615  
 Metal testing, 620  
 Metanephrine test, 621  
 Methemoglobin, 623  
 Methylated septin 9 DNA, 625  
 Metyrapone test, 16  
 Microglobulin, 628  
 Monocytes, 974  
 Mononucleosis rapid test, 630  
*Mycoplasma pneumoniae* antibodies, 631  
 Myoglobin, 635  
 Natriuretic peptides, 636  
 Neuron-specific enolase, 638  
 Neutrophil antibody screen, 639  
 Neutrophil gelatinase-associated lipocalin, 640  
 Newborn metabolic screening, 642  
 Nicotine and metabolites, 646  
 Noninvasive prenatal testing, 217  
 N-telopeptide, 168  
 N-terminal fragment of pro-brain natriuretic peptide, 636  
 5'-Nucleotidase, 649  
 Osmolality, 654  
 Osteocalcin, 168  
 O'Sullivan test, 465  
 Pancreatic enzymes, 662  
 Parathyroid hormone, 675  
 Partial thromboplastin time, activated, 681  
 Parvovirus B19 antibody, 684  
 Pepsinogen, 688  
 Phenylketonuria test, 642  
 Pheochromocytoma suppression and provocative testing, 692  
 Phosphate, 694  
 Phospholipase A2 receptor, 696  
 Phosphorus, 694  
 PI-linked antigen, 697  
 Placental growth factor, 698  
 Plasminogen, 699  
 Plasminogen activator inhibitor I, 700  
 Platelet Aggregation test, 702  
 Antibody detection, 704

## 1004 appendix B: list of tests by type

- Closure time, 709
- Count, 706
- Function assay, 709
- Volume, mean, 712
- Potassium, 724
- Prealbumin, 728
- Pregnancy-associated plasma protein-A, 730
- Pregnancy tests, 515
- Progesterone assay, 733
- Prolactin levels, 737
- Prostate-specific antigen, 743
- Protein, 746
- Protein C, 751
- Protein S, 751
- Prothrombin
  - Fragment, 877
  - Time, 753
- PSA velocity, 743
- QuantiFERON-TB Gold test, 919
- Rabies-neutralizing antibody test, 769
- Red blood cell
  - Count, 770
  - Indices, 772
- Renin assay
  - Plasma, 783
  - Renal vein, 784
- Reticulocyte count, 788
- Rheumatoid factor, 790
- Ribosome P antibodies, 792
- Rubella antibody test, 793
- Rubeola antibody, 615
- SARS, 797
- Septin 9 DNA methylation, 625
- Serotonin, 807
- Sexual assault testing, 809
- Sexually transmitted disease testing, 813
- Sickle cell screen, 820
- Sodium, 835
- Somatomedin C, 541
- Squamous cell carcinoma antigen, 196
- Streptococcus* serologic testing, 850
- Substance abuse testing, 852
- Syphilis detection test, 859
- Testosterone, 861
- Therapeutic drug monitoring, 332
- Thromboelastography, 874
- Thrombosis indicators, 877
- Thyroglobulin, 879
- Thyroid-binding Globulin, 897
  - Inhibitory immunoglobulins, 889
- Thyroid-stimulating Hormone, 886
  - Immunoglobulins, 889
  - Stimulation test, 888
- Thyrotropin receptor antibody, 889
- Thyrotropin-releasing hormone stimulation test, 892
- Thyroxine ( $T_4$ )
  - Free, 894
  - Total, 894
- Thyroxine-binding globulin, 897
- Tissue transglutaminase antibodies, 458
- Toxicology screening, 852
- Toxoplasmosis antibody, 902
- Transferrin, 546
- Transferrin receptor assay, 906
- TRH stimulation test, 886
- Triglycerides, 908
- Triiodothyronine ( $T_3$ ), 910
  - Reverse  $T_3$ , 910
- Troponins, 912
- Tuberculosis testing, 919
- Tumor markers, 194
- Urea nitrogen blood test, 155



- Uric acid, 927
- Uroporphyrinogen-1-synthase, 953
- Viral cultures, 965
- Virus testing, 967
- Vitamin B<sub>12</sub>, 969
- Vitamin D, 971
- Warfarin metabolism genetic test panel, 755
- West Nile virus, 967
- White blood cell count and differential count, 974
- D-Xylose absorption test, 983
- Zinc protoporphyrin, 987

### ELECTRODIAGNOSTIC TESTS

- Caloric study, 192
- Cardiac exercise stress testing, 210
- Contraction stress (fetal), 409
- Electrocardiography, 342
- Electrocorticography, 347
- Electroencephalography, 347
- Electromyography, 350
  - of the pelvic floor sphincter, 353
- Electroneurography, 355
- Electronystagmography, 357
- Electrophysiologic study, 359
- Evoked potential studies, 393
- Fetal
  - Contraction stress test, 409
  - Nonstress test, 415
- Holter monitoring, 510
- Magnetoencephalography, 347
- Microvolt T-wave alternans, 344
- Nonstress (fetal), 415
- Pelvic floor sphincter, electromyography, 353
- Signal-averaged electrocardiography, 344
- Sleep studies, 826

### ENDOSCOPY

- Arthroscopy, 122
- Bronchoscopy, 183
- Capsule endoscopy, 382
- Colonoscopy, 258
- Colposcopy, 262
- Cystoscopy, 309
- Ductoscopy, 336
- Endoscopic retrograde cholangiopancreatography, 364
- Esophagogastroduodenoscopy, 382
- Fetoscopy, 419
- Gastrosocopy, 382
- Hysteroscopy, 533
- Laryngoscopy, 154
- Mediastinoscopy, 618
- Sigmoidoscopy, 258
- Thoracoscopy, 870
- Transesophageal echocardiography, 903

### FLUID ANALYSIS

- Alpha defensin, 39
- Amniocentesis, 49
- Amyloid beta protein precursor, soluble, 58
- Antispermatozoal antibody, 93
- Arthrocentesis with synovial fluid analysis, 118

## 1006 appendix B: list of tests by type

- Beta<sub>2</sub>-microglobulin, 628
- Breast ductal lavage, 179
- Cerebrospinal fluid examination, 576
- Fetal fibronectin, 412
- Genetic testing, 445
- HIV oral testing, 508
- Human papillomavirus, 520
- Hyaluronan binding assay, 803
- Lumbar puncture, 576
- Pancreatic enzymes, 662
- Paracentesis, 671
- Pericardiocentesis, 689
- Placental alpha microglobulin, 52
- SARS, 797
- Semen analysis, 801
- Sexual assault testing, 809
- Sims-Huhner test, 822
- Sperm
  - Chromatin structure assay, 802
  - DNA fragmentation assay, 802
  - DNA integrity, 802
  - Penetration assay, 802
- Substance abuse testing, 852
- Sweat electrolytes test, 857
- Tau protein, 58
- Thoracentesis and pleural fluid analysis, 864
- West Nile virus, 967

### MANOMETRIC TESTS

- Cystometry, 949
- Electrophysiologic study, 359
- Esophageal function studies, 377
- Plethysmography, arterial, 713
- Urethral pressure profile, 949
- Urodynamic studies, 949

### MICROSCOPIC EXAMINATIONS

- Anti-glomerular basement membrane antibodies, 79
- Bioterrorism infectious agents, 141
- Bone marrow biopsy, 160
- Breast cancer
  - Genomics, 454
  - Tumor analysis, 176
- Cathepsin D, 176
- Cervical biopsy, 224
- Chlamydia*, 230
- Colon cancer tumor analysis, 256
- Cutaneous immunofluorescence biopsy, 824
- DNA ploidy status, 176
- Endometrial biopsy, 362
- Estrogen receptor assay, 389
- Gonorrhea culture, 813
- Helicobacter pylori* antibody, 482
- HER 2 protein, 176
- Herpes genitalis, 498
- Ki67 protein, 176
- Liver biopsy, 571
- Lung biopsy, 584
- Lung cancer molecular testing, 588
- Nose culture, 872
- p53 protein, 176
- Pancreatobiliary FISH testing, 665
- Papanicolaou smear, 667
- Parkinson disease testing, 679
- Pleural biopsy, 715
- Progesterone receptor assay, 735

- Renal biopsy, 776
- Sexually transmitted disease cultures, 813
- Skin biopsy, 824
- S-phase fraction, 176
- Streptococcus* screen, 850
- Throat and nose cultures, 872
- Thyroid cancer genomic testing, 455
- Thyroid fine needle aspiration biopsy, 881
- Tuberculosis culture, 917
- Urine culture and sensitivity, 947
- Varicella virus, 958
- Viral cultures, 965
- Wound culture and sensitivity, 981
- Zika virus, 985

### NUCLEAR SCANS

- Bone scan, 165
- Brain scan, 174
- Cardiac nuclear scanning, 207
- Cardiac stress testing, 210
- Gallbladder nuclear scanning, 431
- Gallium scan, 433
- Gastric emptying scan, 437
- Gastroesophageal reflux scan, 441
- Gastrointestinal bleeding scan, 443
- Isonitrile scan, 207
- Liver/spleen scanning, 574
- Lung scan, 590
- Meckel diverticulum nuclear scan, 616
- MUGA scan, 207
- Octreotide scan, 652
- Parathyroid scan, 677
- Positron emission tomography, 719
- ProstaScint scan, 739
- Renal scanning, 779
- Salivary gland nuclear imaging, 795
- Sentinel lymph node biopsy, 805
- Thallium scan, 207
- Thyroid scanning, 883
- Total blood volume, 900
- White blood cell scan, 979

### OTHER STUDIES

- Allergy skin testing, 33
- Bioterrorism infectious agent testing, 141
- Body plethysmography, 762
- Cardiac MRI, 600
- Cell culture drug resistance testing, 216
- Chorionic villus sampling, 241
- Colposcopy, 262
- Ethanol, 391
- Fetal
  - Nonstress test, 415
  - Oxygen saturation monitoring, 417
- Fluorescein angiography, 424
- Gas dilution test, 762
- Genetic testing, 445
- Genomic testing, 453
- Helicobacter pylori* antibody, 482

## 1008 appendix B: list of tests by type

Lactose breath test, 557	Small intestinal bacterial overgrowth, 833
Magnetic resonance	Substance abuse testing, 852
Angiography, 600	Tilt-table test, 359
Imaging, 599	TORCH test, 899
Spectroscopy, 599	Tuberculin test, 915
Venography, 601	Urea breath test, 925
Metal testing, 620	Urine flow studies, 949
Oximetry, 658	Urodynamic studies, 949
Pulmonary function tests, 759	Viral cultures, 965
Sleep studies, 826	

### SPUTUM TESTS

Acid-fast bacillus smear, 917	Culture and sensitivity, 841
Bioterrorism infectious agents, 141	Cytology, 843

### STOOL TESTS

Apt test, 105	<i>Helicobacter pylori</i> antibody, 482
Bioterrorism infectious agents, 141	Lactoferrin, 555
Clostridial toxin assay, 246	Occult blood, 844
Culture, 847	Ova and parasites, 847
Fecal calprotectin, 400	Pancreatic elastase, 660
Fecal fat, 402	Stool for leukocytes, 849
	Viral cultures, 965

### ULTRASOUND TESTS

Abdominal ultrasound, 1	Obstetric ultrasonography, 685
Amniotic fluid analysis, 49	Pelvic ultrasonography, 685
Breast sonogram, 181	Prostate/rectal sonogram, 741
Carotid artery duplex scanning, 214	Scrotal ultrasound, 799
Echocardiography, 339	Testicular ultrasound, 799
Endourethral urologic ultrasound, 1	Thyroid ultrasound, 891
Fetal	Transesophageal echocardiography, 903
Biophysical profile, 406	Transrectal ultrasonography, 741
Nuchal translucency, 685	Transthoracic echocardiography, 339
Intravascular ultrasound, 543	Vascular ultrasound studies, 960
IUD localization, 686	
Kidney sonogram, 1	
Liver and pancreaticobiliary system ultrasonography, 1	

## URINE TESTS

- Adrenocorticotrophic hormone stimulation test with metyrapone, 16  
 Aldosterone, 25  
 Amino acid profiles, 45  
 Amylase, 56  
 Appearance and color, 931  
 Bence Jones protein, 134  
 Beta<sub>2</sub>-microglobulin, 628  
 11 Beta-prostaglandin F(2) alpha, 136  
 Bioterrorism infectious agents, 141  
 Bladder cancer markers, 145  
 Bone turnover markers, 168  
 Calcium, 189  
 Cortisol, 290  
 Creatinine clearance, 303  
 Crystals, 936  
 Culture and sensitivity, 947  
 Delta-aminolevulinic acid, 321  
 Dexamethasone suppression test, 324  
 Epithelial casts, 937  
 Estimated glomerular filtration rate, 303  
 Estriol excretion, 386  
 Estrogen fractions, 386  
 Ethanol, 391  
 Fatty casts, 936  
 Flow studies, 949  
 Fluorescence *in situ* hybridization, 145  
 Fractional excretion of sodium, 837  
 Glucose, 934  
 Glucose tolerance test, 467  
 Granular casts, 936  
 HIV urine testing, 508  
 Hyaline casts, 936  
 17-Hydroxycorticosteroids, 526  
 5-Hydroxyindoleacetic acid, 528  
 Ketones, 935  
 17-Ketosteroids, 862  
 Leucine aminopeptidase, 561  
 Leukocyte esterase, 935  
 Liquid biopsy, 569  
 Metyrapone test, 16  
 Microalbumin, 626  
 Nicotine and metabolites, 646  
 N-telopeptide, 168  
 Nitrites, 935  
 Odor, 931  
 Osmolality, 656  
 Pepsinogen, 688  
 pH, 931  
 Phenylketonuria test, 642  
 Porphyrins and porphobilinogens, 717  
 Potassium, 724  
 Prealbumin, 728  
 Pregnancy tests, 515  
 Pregnanediol, 731  
 Protein, 746  
 Pyridinium, 168  
 Red blood cells and casts, 937  
 Sodium, 835  
 Specific gravity, 934  
 Substance abuse testing, 852  
 Toxicology screening, 852  
 Tubular casts, 937  
 Uric acid, 927  
 Urinalysis, 930  
 Urinary methylmalonic acid, 969  
 Urinary stone analysis, 945

## 1010 appendix B: list of tests by type

- Urine culture and sensitivity, 947
- Vanillylmandelic acid and catecholamines, 954
- Viral cultures, 965
- Waxy casts, 936
- White blood cells and casts, 937
- D-Xylose absorption test, 983

### X-RAY EXAMINATIONS

- Arteriography, 114
- Barium
  - Enema, 128
  - Swallow, 131
- Bone
  - Densitometry, 157
  - X-ray, 172
- Cardiac catheterization, 201
- Chest x-ray, 227
- Computed tomography
  - Arteriography, 273
  - Colonography, 273
  - of the abdomen, 272
  - of the adrenals, 272
  - of the brain, 278
  - of the chest, 281
  - of the heart, 284
  - of the kidney, 272
- Cystography, 307
- Dental x-rays, 322
- Digital subtraction angiography, 114
- Hysterosalpingography, 531
- Intravenous pyelography, 764
- Kidney, ureter, bladder x-ray study, 650
- Magnetic resonance imaging, 599
- Mammography, 606
- Myelography, 632
- Obstruction series, 650
- Percutaneous transhepatic cholangiography, 364
- Positron emission tomography, 719
- Pulmonary angiography, 757
- Pyelography, 764
- Renal angiography, 114
- Retrograde pyelography, 764
- Sialography, 818
- Skull x-ray, 825
- Small bowel follow-through, 830
- Spinal x-rays, 839
- Swallowing examination, 855
- Upper gastrointestinal x-ray study, 922
- Venography of lower extremities, 963
- Videofluoroscopy, 855
- Virtual colonoscopy, 273

## appendix C: disease and organ panels

This appendix includes groupings of tests commonly used either for screening or to evaluate disease situations. These panels may be modified or expanded in different clinical settings.

### Anemia

- Complete blood count (CBC), 269
- Macrocytic anemia
  - Folate, 426
  - TSH, 886
  - Vitamin B<sub>12</sub>, 969
- Microcytic anemia
  - Reticulocyte index, 788
  - Iron panel, 546
- Normocytic anemia
  - Reticulocyte index, 788
  - Hemolysis profile, 1012
- Red blood cell (RBC) indices, 772
- Reticulocyte count, 788

### Arthritis

- Antinuclear antibody (ANA), 86
- C-reactive protein, 295
- ESR (sedimentation rate), 373
- Rheumatoid factor, 790
- Uric acid, 927

### Basic metabolic panel (BMP)

- Blood urea nitrogen (BUN), 155
- Calcium, 189
- Carbon dioxide, 197
- Chloride, 233
- Creatinine, 301
- Glucose, 462
- Potassium, 724
- Sodium, 835

### Bone/joint

- Albumin, 746
- Alkaline phosphatase, 29
- Calcium, 189
- Osteocalcin, 168
- Phosphorus, 694
- Protein, total, 746
- Uric acid, 927

### Cardiac injury

- Creatine kinase (CK), 297
- CK-MB, 297
- Myoglobin, 635
- Troponin I, 912

### Coagulation screening

- Partial thromboplastin time (PTT), 681
- Platelet count, 706
- Prothrombin time (PT), 753

### Coma

- Alcohol, 391
- Ammonia, 47
- Anion gap, 63
- Arterial blood gases, 106
- Basic metabolic panel (BMP), 1011
- Calcium (total and ionized), 189
- Ethyl alcohol, 391
- Lactic acid, 551
- Osmolality (serum), 654
- Salicylate, 333
- Toxicology screen, 852

### **Comprehensive metabolic panel (CMP)**

Albumin, 746  
Alkaline phosphatase, 29  
Aspartate aminotransferase (AST, SGOT), 125  
Bilirubin, 137  
Blood urea nitrogen (BUN), 155  
Calcium, 189  
Carbon dioxide, 197  
Chloride, 233  
Creatinine, 301  
Glucose, 462  
Potassium, 724  
Protein, total, 746  
Sodium, 835

### **Diabetes mellitus management**

Anion gap, 63  
Basic metabolic panel (BMP), 1011  
Hemoglobin A<sub>1c</sub>, 471  
Lipid profile, 565

### **DIC**

Complete blood count (CBC), 269  
Fibrin split products, 319  
Fibrinogen, 422  
Partial thromboplastin time (PTT), 681  
Platelet count, 706  
Prothrombin time (PT), 753

### **Electrolyte panel**

Carbon dioxide, 197  
Chloride, 233  
Potassium, 724  
Sodium, 835

### **General health**

Complete blood count (CBC), 269

Comprehensive metabolic panel (CMP), 270  
Gamma-glutamyl transferase (GGT), 435  
Lactic dehydrogenase (LDH), 553  
Lipid profile, 565  
Thyroid-stimulating hormone (TSH), 886  
Uric acid, 927

### **Hemolysis profile**

Antiglobulin, 287, 289  
Bilirubin, 137  
Complete blood count (CBC), 269  
Haptoglobin, 479  
Hemoglobin, 488  
Lactic dehydrogenase (LDH), 553  
Reticulocyte count, 788

### **Hepatic function panel**

Alanine aminotransferase (ALT)/serum glutamic-pyruvic transaminase (SGPT), 21  
Albumin, 746  
Alkaline phosphatase, 29  
Aspartate aminotransferase (AST)/serum glutamic-oxaloacetic transaminase (SGOT), 125  
Bilirubin, direct and total, 137  
Gamma-glutamyl transferase (GGT), 435  
Protein, total, 746  
Prothrombin time, 753

### **Hepatitis, acute**

Hepatitis A antibody IgM, 493  
Hepatitis B, 493  
Core antibody IgM, 493



Surface antigen, 493  
Hepatitis C antibody, 493

## HIV

CD4 and CD8, 219  
Complete blood count (CBC), 269  
HIV antibody with Western blot confirmation, 506

## Hypertension

Basic metabolic panel (BMP), 1011  
Cortisol, urinary free, 290  
Metanephrines, urinary, 954  
Renin, 783  
Thyroid screening panel, 1014  
Urinalysis, 930

## Lipid panel

HDL cholesterol, 565  
Total cholesterol, 235  
Triglycerides, 908

**Liver function panel** – See  
Hepatic function panel, 1012

## Obstetric

Blood  
Screen, 152  
Type, 152  
Complete blood count (CBC), 269  
Hepatitis B surface antigen, 494  
Rubella antibody, 793  
Syphilis test (RPR, VDRL), 859

## Pancreatic

Amylase, 56  
Calcium (total and ionized), 189  
Glucose, 462  
Lipase, 562  
Triglycerides, 908

## Parathyroid

Albumin, 746  
Alkaline phosphatase, 29  
Calcium (total and ionized), 189  
Calcium, urinary, 189  
Creatinine, 301  
Magnesium, 597  
Parathyroid hormone (PTH), 677  
Phosphorus, 694  
Protein, total, 746  
PTH, 677

## Prenatal

ABO and Rh typing, 152  
Antibody screen, 289  
Blood urea nitrogen (BUN), 155  
Cervical cultures for sexually transmitted diseases, 231  
Complete blood count (CBC), 269  
Creatinine, 301  
Cytomegalovirus (CMV), 317  
Glucose, 462  
Hepatitis B surface antigen, 494  
Herpes simplex I and II, 498  
Pap smear, 667  
Rubella titer, 793  
Thyroxine, free  $T_4$ , 894  
Toxoplasmosis antibody, 902  
Uric acid, 927  
Urinalysis, 930  
Urine culture, 947  
VDRL, 859

## Renal panel

Albumin, 746  
Blood urea nitrogen (BUN), 155  
Calcium, 189

## 1014 appendix C: disease and organ panels

Carbon dioxide, 197  
Chloride, 233  
Complete blood count  
(CBC), 269  
Creatinine  
  Blood, 301  
  24-hour urine, 303  
Glucose, 934  
Magnesium, 597  
Phosphorus, 694  
Potassium, 724, 727  
Sodium, 835

### Thyroid screening panel

Thyroid-stimulating hormone  
(TSH), 886, 888  
Thyroxine (free  $T_4$ ), 894

### TORCH antibody panel

Cytomegalovirus antibody, 317  
Herpes simplex antibody, 498  
Rubella antibody, 793  
Toxoplasmosis antibody, 902

### Toxicology screening (urine)

Amphetamines, 852  
Barbiturates, 852  
Benzodiazepines, 852  
Cocaine metabolites, 852  
Marijuana metabolites, 852  
Methadone, 853  
Opiate metabolites, 852  
Phencyclidine, 852  
Propoxyphene, 852

## appendix D: symbols and units of measurement

<	less than
≤	less than or equal to
>	greater than
≥	greater than or equal to
c (prefix)	centi ( $10^{-2}$ )
C	celsius
cc	cubic centimeter
cg	centigram
cm	centimeter
cm H <sub>2</sub> O	centimeter of water
cu	cubic
d (prefix)	deci ( $10^{-1}$ )
dL	deciliter (100 mL)
f (prefix)	femto ( $10^{-15}$ )
fL	femtoliter
fmol	femtomole
g	gram
hr	hour
IU	international unit
ImU	international milliunit
I $\mu$ U	international microunit
k (prefix)	kilo ( $10^3$ )
kat	katal
kg	kilogram (1000 grams)
L	liter
m (prefix)	milli ( $10^{-3}$ )
m	meter
m <sup>2</sup>	square meter
m <sup>3</sup>	cubic meter
mcg	microgram
mEq	milliequivalent
mEq/L	milliequivalent per liter
mg	milligram (1/1000 gram)
min	minute
mL	milliliter
mm	millimeter (1/10 centimeter)
mm <sup>3</sup>	cubic millimeter
mM	millimole
mm Hg	millimeter of mercury
mm H <sub>2</sub> O	millimeter of water

**1016** appendix D: symbols and units of measurement

mmol	millimole
mol	mole
mOsm	milliosmole
m $\mu$	millimicron
mU	milliunit
mV	millivolt
$\mu$ (prefix)	micro ( $10^{-6}$ )
$\mu^3$	cubic micron
$\mu$ kat	microkatal
$\mu$ L	microliter
$\mu$ m	micrometer
$\mu$ m <sup>3</sup>	cubic micrometer
$\mu$ mol	micromole
$\mu$ U	microunit
n (prefix)	nano ( $10^{-9}$ )
ng	nanogram
nkat	nanokatal
nm	nanometer
nmol	nanomole
p (prefix)	pico ( $10^{-12}$ )
Pa	pascal
pg	picogram
pL	picoliter
pm	picometer
pmol	picomole
sec	second
SI units	International System of Units
U	unit
yr	year

## bibliography

- Ashoor G, Syngelaki A, Nicolaides KH, et al. Trisomy 13 detection in the first trimester of pregnancy using a chromosome-selective cell-free DNA analysis method. *Ultrasound Obstet Gynecol* 41(1):21–25, 2013.
- Beta J, Bredaki FE, Calvo JR, et al. Maternal serum  $\alpha$ -fetoprotein at 11–13 weeks' gestation in spontaneous early preterm delivery. *Fetal Diagn Ther* 30:88–93, 2011.
- Brenner H, Chang-Claude J, Seiler CM, et al. Protection from colorectal cancer after colonoscopy. *Ann Intern Med* 154:22–30, 2011.
- Centers for Disease Control and Prevention. Seasonal influenza. <http://www.cdc.gov/flu>.
- Chang H, Shin BK, Kim A, et al. DNA methylation analysis for the diagnosis of thyroid nodules—a pilot study with reference to BRAF<sup>V600E</sup> mutation and cytopathology results. *Cytopathology* 27:122–130, 2016.
- Chapman PB, Hauschild A, Robert C, et al. Improved survival with vemurafenib in melanoma with BRAF V600E mutation. *N Engl J Med* 364(26):2507–2516, 2011.
- Cheng JM, Suoniemi M, Kardys I, et al. Plasma concentrations of molecular lipid species in relation to coronary plaque characteristics and cardiovascular outcome: results of the ATHEROREMO-IVUS study. *Atherosclerosis* 243:560–566, 2015.
- Costa CR, Johnson AJ, Naziri Q, et al. Efficacy of erythrocyte sedimentation rate and C-reactive protein level in determining periprosthetic hip infections. *Am J Orthop (Belle Mead NJ)* 41(4):160–165, 2012.
- Felker GM, Fiuzat M, Shaw LK, et al. Galectin-3 in ambulatory patients with heart failure: Results from the HF-ACTION study. *Circulation* 125(1):72–78, 2012.
- Kägi G, Bhatia KP, Tolosa E. The role of DAT-SPECT in movement disorders. *J Neurol Neurosurg Psychiatry* 81(1):5–12, 2010.
- Kitchener HC, Gilham C, Sargent A, et al. A comparison of HPV DNA testing and liquid based cytology over three rounds of primary cervical screening: extended follow-up in the ARTISTIC trial. *Eur J Cancer* 47:864–871, 2011.
- Klein EA, Cooperberg MR, Magi-Galuzzi C, et al. A 17-gene assay to predict prostate cancer aggressiveness in the context of Gleason grade heterogeneity, tumor multifocality, and biopsy undersampling. *Eur Urol* 66(1):550–560, 2014.
- Kurtz SM, Lau E, Watson H, Schmier JK, Parvizi J. Economic burden of periprosthetic joint infection in the United States. *J Arthroplasty* 27:61–65, 2012.
- Lea DH, Skirton H, Read CY, et al. Implications for educating the next generation of nurses on genetics and genomics in the 21st century. *J Nurs Schol* 43(1):3–12, 2011.
- Li N, Franceschi S, Howell-Jones R, et al. Human papillomavirus type distribution in 30,848 invasive cervical cancers worldwide: variation by geographical region, histological type and year of publication. *Int J Cancer* 128:927–935, 2011.

## 1018 bibliography

- Maiz N, Wright D, Ferreira AF, et al. A mixture model of ductus venosus pulsatility index in screening for aneuploidies at 11-13 weeks' gestation. *Fetal Diagn Ther* 31:221-229, 2012.
- Nicolaides KH, Syngelaki A, Ashoor G, et al. Noninvasive prenatal testing for fetal trisomies in a routinely screened first-trimester population. *Am J Obstet Gynecol* 207:374-376, 2012.
- Nicolaides KH, Wright D, Poon LC. First-trimester contingent screening for trisomy 21 by biomarkers and maternal blood cell-free DNA testing. *Ultrasound Obstet Gynecol* 42:41-50, 2013.
- Norton M, Brar H, Weiss J, et al. Non-Invasive Chromosomal Evaluation (NICE) Study: results of a multicenter, prospective, cohort study for detection of fetal trisomy 21 and trisomy 18. *Am J Obstet Gynecol* 7(2):137-138, 2012.
- Oduyebo T, Igbinsola I, Petersen EE, et al. Update: interim guidance for health care providers caring for pregnant women with possible Zika virus exposure—United States, July 2017. *MMWR Morb Mortal Wkly Rep* 66(29):781-793, 2017.
- Pagana KD, Pagana TJ. *Mosby's manual of diagnostic and laboratory tests*, ed 6. Mosby, 2018, St. Louis.
- Pandya P, Wright D, Syngelaki A, et al. Maternal serum placental growth factor in prospective screening for aneuploidies at 8-13 weeks' gestation. *Fetal Diagn Ther* 31:87-93, 2012.
- Quintero E, Castells A, Bujanda L, et al. Colonoscopy versus fecal immunochemical testing in colorectal cancer screening. *N Engl J Med* 366:697-706, 2012.
- Schlumberger W, Hornig N, Lange S, et al. Differential diagnosis of membranous nephropathy with autoantibodies to phospholipase A2 receptor 1. *Autoimmun Rev* 13(2):108-113, 2014.
- Waggoner JJ, Pinsky BA. Zika virus: diagnostics for an emerging pandemic threat. *J Clin Microbiol* 54(4):860-867, 2016.
- Yu J, Pan W, Shi R, et al. Ceramide is upregulated and associated with mortality in patients with chronic heart failure. *Can J Cardiol* 31(3):357-363, 2015.
- Ziai J, Hui P. BRAF mutation testing in clinical practice. *Expert Rev Mol Diagn* 12(2):127-138, 2012.

Note: Page numbers followed by *f* indicate figures, *t* indicate tables and *b* indicate boxes.

## A

- AB index, 962  
 AB<sub>42</sub>/AB<sub>40</sub> ratio, 58  
 Abdomen, computed tomography of, 272–277, 277*t*  
 Abdominal aorta, ultrasound of, 1, 4  
 Abdominal cavity, ultrasound of, 4  
 Abdominal scintigraphy, 443–444, 444*t*  
 Abdominal sonogram, 1–4  
 Abdominal ultrasound, 1–4, 2*f*  
 ABO system, 152, 152*t*  
 Absolute neutrophil count (ANC), 976  
 Abuse testing, substance, 852–854, 853*t*  
 Academy rash, 684  
 Acetaminophen, 333*t*  
 Acetazolamide (Diamox), 174  
 Acetyl salicylic acid (ASA), 709  
 Acetylcholine (ACh), 5  
 Acetylcholine receptor (AChR), 5  
 Acetylcholine receptor antibody panel (AChR Ab), 5–6, 6*t*  
 Acetylcholinesterase, 238–240, 240*t*  
 AChR-binding antibody, 5  
 AChR-blocking antibody, 5  
 AChR-modulating antibody, 5  
 Acid clearing, 378  
 Acid perfusion, 378  
 Acid phosphatase, 7–8, 8*t*  
 Acid reflux, with pH probe, 377–378  
 Acid-base balance, 724  
   in K excretion, 727  
 Acid-base disturbances, 108–109*t*  
 Acidemia, 931–933  
 Acid-fast bacillus (AFB) smear, 917–918  
 Acidosis  
   hypoxia and, 417  
   lactic, 551, 552*t*  
   metabolic, 108–109*t*  
   respiratory, 108–109*t*  
 aCL antibodies, 65–66, 66*t*  
 Acquired AAT deficiency, 37  
 Acquired AT-III deficiency, 97  
 Acquired deficiency, of alpha<sub>1</sub>-antitrypsin, 37  
 Acquired immune deficiencies, immunoglobulins in, 537  
 Acquired immunodeficiency syndrome (AIDS)  
   cell surface immunophenotyping for, 219  
   serology for, 506–509  
   T-lymphocyte cell markers for, 219–221, 219*t*, 221*t*  
 Actigraphy, 828  
 Activated clotting time (ACT), 9–10, 10*t*  
 Activated coagulation time, 9–10, 10*t*  
 Activated partial thromboplastin time (APTT), 9, 681–683, 683*t*  
 Activated protein C (APC) resistance test, 397  
 Activating procedures, EEG and, 349  
 Active enzyme, pepsinogen and, 688  
 Active pulmonary sarcoidosis, angiotensin-converting enzyme in, 61  
 Active ulcer disease, *Helicobacter pylori* testing in, 482  
 Acute coronary syndromes, PAPP-A in, 730  
 Acute EBV infection, 368  
 Acute extrahepatic obstruction, aspartate aminotransferase in, 125  
 Acute hepatitis, aspartate aminotransferase in, 125  
 Acute infections, in parvovirus B19 antibody, 684  
 Acute intermittent porphyria (AIP), 321, 953  
 Acute iron poisoning, serum iron in, 547  
 Acute kidney injury, neutrophil gelatinase-associated lipocalin for, 640  
 Acute renal failure, neutrophil gelatinase-associated lipocalin for, 640  
 Acute stroke, ischemia-modified albumin in, 550  
 Acute tubular necrosis, 837  
 Acutely ill patients, total blood volume for, 900  
 Acute-phase protein, 373  
 Acute-phase reactant protein  
   ferritin and, 404  
   fibrinogen as, 422  
 Addison disease  
   17-OHCS in, 526  
   adrenocorticotrophic hormone in, 11  
   adrenocorticotrophic hormone stimulation test with cosyntropin in, 13  
   cortisol and, 290  
   21-hydroxylase antibodies in, 530  
   in metyrapone test, 16  
   SIADH and, 72

## 1020 index

- Addisonian crisis, in metyrapone test, 16
- Adenohypophysis, prolactin levels and, 737
- Adenoma, 677
  - aldosterone in, 25
- Adenosine, 211
- ADH stimulation test, 71–72
- ADH suppression test, 72
- Adrenal gland
  - hyperfunctioning of, 17-OCHS and, 526
  - hypofunctioning of, 17-OCHS and, 526
- Adrenal hyperplasia, metyrapone test in, 16
- Adrenal insufficiency
  - adrenocorticotropic hormone stimulation test with cosyntropin in, 13–14
  - metyrapone test in, 16
- Adrenal nodular hyperplasia, bilateral, aldosterone in, 25
- Adrenal steroid precursors, 18–19, 18–19t
- Adrenal tumors, adrenal steroid precursors in, 19
- Adrenaline, WBC levels and, 976
- Adrenocorticotropic hormone (ACTH), 11–12, 12t, 324
  - stimulation test
    - with cosyntropin, 13–15
    - with metyrapone, 16–17
  - suppression test, 324–326, 326t
- Aedes* mosquito, Zika virus in, 985
- Afirma gene expressive classifier (GEC), 456
- Agatston score, 284–285, 285t
- Age-adjusted PSA, 743
- Age-related macular degeneration (ARMD), risk analysis of, 20
- Agglutination, activity of, 702
- Agglutinins, 774
  - febrile, 399, 399t
- AIDS (acquired immunodeficiency syndrome)
  - cell surface immunophenotyping for, 219
  - serology for, 506–509
  - T-lymphocyte cell markers for, 219–221, 219t, 221t
- Air-contrast barium enema, 128
- Air-contrast upper GI study, 923
- Airflow assessment, pulmonary function tests and, 759
- Airflow monitors, 826
- AKTI, 588
  - testing for, 588–589
- Alanine aminotransferase (ALT), 21–22, 22t
- Albumin, 746–750, 750t
  - in cerebrospinal fluid, 578
  - ischemia-modified, 550
  - ratio, 747
- Albumin gradient, 865–866
- Alcohol
  - ethyl, catecholamines and, 956
  - levels, ethanol (EtOH), 391–392
  - testing, 852
  - testosterone and, 862
- Aldolase, 23–24, 24t
- Aldosterone, 25–28, 27–28t, 724, 727, 835
  - Aldosterone stimulation test, 26
  - Aldosterone suppression test, 26
  - Aldosterone/renin ratio, 783–784
- Aldosteronism, aldosterone in
  - primary, 25–28
  - secondary, 25, 27–28
- ALK, 588
  - gene, 195–196t
  - testing for, 588–589
- Alkaline, acid phosphatase and, 7
- Alkaline phosphatase (ALP), 29–30, 30t
  - bone-specific, serum, 169
  - paracentesis and, 671–672
- Alkalosis
  - metabolic, 108–109t
  - respiratory, 108–109t
- Alkalotic states, 727
- Allergy blood testing, 31–32, 32t
- Allergy skin testing, 33–36, 36t
  - in IgE antibodies, 31
  - intradermal method in, 35
  - patch method in, 35
  - prick-puncture method in, 35
- Alpha defensin test, 39–40
- Alpha<sub>1</sub>-antitrypsin (A<sub>1</sub>AT), 37–38, 38t
  - phenotyping, 37–38, 38t
- Alpha-1-microglobulin, 628–629, 629t
- Alpha-fetoprotein (AFP), 41–42, 42t, 195–196t
- Alpha-synuclein, 679
- Aluminum, 43–44, 44t
- Aluminum-based phosphate binder gels, 43
- Alveolar to arterial O<sub>2</sub> difference (A-a gradient), 106, 110
- Alzheimer disease (AD), 720
  - amyloid beta protein precursor in, 58
- American Cancer Society, mammography screening guidelines, 608–609t
- American College of Obstetricians and Gynecologists (ACOG), 217
  - mammography screening guidelines, 608–609t
- American Society of Breast Surgeon, mammography screening guidelines, 608–609t
- Amine precursor uptake and decarboxylation (APUD) cells, 638



- Amino acids, 45  
 profiles, 45–46, 46*t*  
 screen, 45–46, 46*t*
- 40-amino-acid peptide (AB<sub>40</sub>), 58
- 42-amino-acid peptide (AB<sub>42</sub>), 58
- Aminolevulinic acid (ALA), 321, 321*t*, 717
- Aminophylline, 333*t*  
 catecholamines and, 956
- Amino-terminal propeptide of type I procollagen (PINP), 168, 170
- Ammonia, 47–48, 48*t*  
 paracentesis and, 671–672
- Amniocentesis, 49–55, 49*t*, 53*f*, 55*t*
- Amnio-dye test, 51–52
- Amniotic fluid  
 analysis, 49–55, 49*t*, 55*t*  
 for infection, assessment of, 51  
 optical density of, measurement of, 50  
 volume, 407
- Amniotic fluid index (AFI), 407–408
- Amphetamines, 853*t*
- Amylase, 56–57, 57*t*, 662–664, 664*t*  
 paracentesis and, 671–672  
 in pleural fluid analysis, 866
- Amylase/creatinine clearance ratio, 57
- Amyloid, PET scanning with, 720
- Amyloid beta protein precursor, soluble (sBPP, APP), 58–59, 59*t*
- Anabasine, 647
- Anabolic steroids  
 thyroxine-binding globulin and, 897  
 total thyroxine (T<sub>4</sub>) and, 895  
 triiodothyronine (T<sub>3</sub> radioimmunoassay) and, 911
- Anaerobic glycolysis, fetal hypoxia causing, 417
- Anaerobic metabolism, of glucose, 551
- Anal canal culture, 815
- Anaphylactic reaction, in allergen, 31
- ANAs (antinuclear antibodies), 67, 74–75, 86–89, 86–87*b*, 87*t*, 88*f*, 89*t*, 95, 867
- Anatomic abnormalities, in fetal lung maturity test, 51
- Androgens  
 testosterone and, 861–862  
 thyroxine-binding globulin and, 897  
 total thyroxine (T<sub>4</sub>) and, 895  
 triiodothyronine (T<sub>3</sub> radioimmunoassay) and, 911
- Androstenediones (AD), 18–19, 18–19*t*
- Anemia  
 categorization of, 773*t*  
 2,3-diphosphoglycerate and, 328  
 erythropoietin for, 375  
 hematocrit in, 485
- Anemia (*Continued*)  
 hemoglobin in, 488, 774  
 iron deficiency, 404, 547, 987  
 pernicious, 545, 688, 969  
 red blood cell count for, 771*t*  
 red blood cell distribution width for, 774  
 red blood cell indices for, 772  
 reticulocyte count for, 788  
 sickle cell, 643, 820  
 tests for, 1011–1014  
 total blood volume for, 900  
 transferrin receptor (TfR) assay for, 906
- Angiography. *See* Arteriography (angiography)
- Angiotensin, 60, 60*t*
- Angiotensin-converting enzyme (ACE), 60–62, 61–62*t*
- Anion gap (AG), 63–64, 64*t*
- Anisocytosis, 149, 774
- Ankle-brachial index (ABI), 713–714, 962
- Anoscopy, 258–261
- Antegrade pyelography, 764–768
- Anterior pituitary gland  
 adenohypophysis, prolactin levels and, 737  
 adrenocorticotropic hormone in, 11
- Anthrax, bioterrorism infectious agents testing for, 141, 142*t*
- Anti-AChR antibody, 5–6, 6*t*
- Antiarrhythmic drugs, effectiveness of, 359
- Antibiotics, 119  
 blood cultures and, 147  
 thromboelastography and, 874–875  
 WBC levels and, 976
- Antibody screening, blood, 289
- Antibody titer, Epstein-Barr virus, 368–370, 369*t*
- Antibody(ies)  
 anticardiolipin, 65–66, 66*t*  
 anticomplement, 67  
 antichitobioside carbohydrate, 80  
 antichromatin, 68–69, 69*t*  
 anticyclic citrullinated peptide, 70, 70*t*  
 antifungal, 428–429, 429*t*  
 antiglomerular basement membrane, 79  
 antiglycan, 80, 80*t*  
 anti-*Helicobacter pylori*, 482–484  
 antihistidyl transfer synthase, 75  
 antihistone, 68–69, 69*t*  
 antilaminaribioside carbohydrate, 80  
 antiliver/kidney microsomal type 1, 81, 81*t*  
 anti-LKM-1, 81, 81*t*  
 antimannobioside carbohydrate, 80  
 antimitochondrial, 82, 82*t*  
 antimyocardial, 83, 83*t*  
 antineutrophil cytoplasmic, 84–85, 84–85*t*

## 1022 index

- Antibody(ies) (*Continued*)
- antinuclear, 67, 74–75, 86–89, 86–87*b*, 87*t*, 88*f*, 89*t*, 95, 867
  - antinucleosome, 68–69, 69*t*
  - antiparietal cell, 90, 90*t*
  - antiphospholipid, 65–66, 66*t*
  - anti-*Saccharomyces cerevisiae*, 80
  - antiscleroderma, 91
  - antismooth muscle, 92, 92*t*
  - antispermatozoal, 93–94
  - anti-SS-A, 95–96
  - anti-SS-B, 95–96
  - anti-SS-C, 95–96
  - antithyroglobulin, 99, 99*t*
  - antithyroid peroxidase, 100, 100*t*
  - antivinculin, 316
  - in blood group, 152, 152*t*
  - to double-stranded DNA, 74, 74*t*
  - endomyxial, 458–459, 458*t*
  - Epstein-Barr virus, 368–370, 369*t*
    - for mononucleosis, 368
  - to extractable nuclear antigens, 75–76, 75–76*t*
  - febrile, 399, 399*t*
  - gliadin, 458–459, 458*t*
  - glutamic acid decarboxylase, 327, 327*t*
  - granulocyte, 639
  - histone, 68
  - human T-cell lymphotropic virus (HTLV) I/II, 525
  - 21-hydroxylase, 530
  - intrinsic factor, 545, 545*t*
  - islet cell, 327
  - neutrophil, 639
  - phospholipid, 65
  - polymorphonucleocyte, 639
  - RNA polymerase III, 91
  - scleroderma, 91
  - Sjögren, 95–96
  - skin biopsy, 824
  - sperm, 93–94
  - test
    - for HIV, 506–509
    - rabies-neutralizing, 769
    - rubella, 793–794, 793*t*
    - for varicella zoster virus, 958–959
  - thyroid-stimulating hormone, 889–890
  - thyrotropin receptor, 889–890
  - tissue transglutaminase, 458–459, 458*t*
  - toxoplasmosis, 902
  - tuberculosis, 919–921
- Anticardiolipin antibodies (ACA), 65–66, 66*t*
- Anticentromere antibody test, 67
- Antichitobioside carbohydrate antibody (ACCA), 80
- Antichromatin antibody test, 68–69, 69*t*
- Anticoagulation, INR according to indication for, 754*t*
- Anticonvulsants
- testosterone and, 862
  - WBC levels and, 976
- Anticyclic citrullinated peptide antibody, 70, 70*t*
- Anticytoplasmic antibody, 92
- Antideoxyribonuclease-B titer (anti-DNase-B, ADNase-B, ADB), 850–851, 851*t*
- Anti-deoxyribonucleic acid antibodies, 74, 74*t*
- Antidiuretic hormone (ADH), 71–73, 73*t*, 835
- Anti-DNA antibody test, 74, 74*t*
- Anti-double-stranded DNA, 74, 74*t*
- Anti-ds-DNA, 74, 74*t*
- Anti-EA (EA-D or EA-R) antibodies, 369
- Antixtractable nuclear antigens (anti-ENAs), 75–76, 75–76*t*, 95
- Antifactor Xa (anti-Xa), 77–78
- Antifungal antibodies, 428–429, 429*t*
- Anti-GBM antibody, 79
- Antigens
- assay, 97–98, 98*t*
  - in blood group, 152, 152*t*
  - platelet-specific, 704
  - Rh, red blood cell and, 153
  - Streptococcus* group B, 851
  - white blood cell, 518–519
- Antigliadin, 458
- Antiglomerular basement membrane antibodies (AGBM), 79
- Antiglycan antibodies, 80, 80*t*
- Antigranulocyte antibodies, 639
- Anti-HCV antibodies, 496
- Anti-*Helicobacter pylori* antibody, 482–484
- Antihemophilic factor, 249–250*t*, 252*t*
- Antihistamines, WBC levels and, 976
- Antihistidyl transfer synthase (anti-Jo-1), 75–76, 75–76*t*
- Antihistone antibody test (anti-HST, AHA), 68–69, 69*t*
- Anti-La antibodies, 95–96
- Antilaminaribioside carbohydrate antibody (ALCA), 80
- Antiliver/kidney microsomal type 1 antibodies, 81, 81*t*
- Anti-LKM-1 antibodies, 81, 81*t*
- Antimannobioside carbohydrate antibody (AMCA), 80
- Antimetabolites, WBC levels and, 976
- Antimitochondrial antibody (AMA), 82, 82*t*
- Anti-mutated citrullinated vimentin (anti-MCV), 70, 70*t*
- Antimyocardial antibody (AMA), 83, 83*t*
- Antineutrophil antibodies, 639
- Antineutrophil cytoplasmic antibody (ANCA), 84–85, 84–85*t*

- Antinuclear antibodies (ANAs), 67, 74–75, 86–89, 86–87*b*, 87*t*, 88*f*, 89*t*, 95, 86*t*
- Antinuclear ribonucleoprotein (anti-RNP) antibody, 75
- Antinuclear Smith (anti-SM) antibody, 75
- Antinucleosome antibody (anti-NCS), 68–69, 69*t*
- Antiparietal cell antibody (APCA), 90, 90*t*
- Antiphospholipid antibodies, 65–66, 66*t*  
activated partial thromboplastin time and, 68*t*
- Antiphospholipid syndrome (APS),  
antcardiolipin antibodies in, 65
- Antiplatelet antibody detection, 704–705, 705*t*
- Antiplatelet drugs, thromboelastography and, 874–875
- Antiretroviral therapy, based on viral load and CD4 count, 504*t*
- Antiribonucleoprotein (anti-RNP), 75–76, 75–76*t*
- Anti-ribosome P antibodies, 792, 792*t*
- Anti-Ro antibodies, 95–96
- Antirubella antibody testing, 794
- Anti-*Saccharomyces cerevisiae* antibody (ASCA), 80
- Antiscleroderma antibody, 91
- Anti-Smith (anti-SM), 75–76, 75–76*t*
- Antismooth muscle antibody (ASMA), 92, 92*t*
- Antisperm antibodies, 93–94
- Antispermatozoal antibody, 93–94
- Anti-SS-A (RO), 95–96
- Anti-SS-B (LA), 95–96
- Anti-SS-C antibodies, 95–96
- Antistreptolysin O (ASO) antibody, 850–851
- Antistreptolysin O (ASO) titer, 850–851, 851*t*
- Antistriated muscle antibody, 5
- Antithrombin activity, 97–98, 98*t*
- Antithrombin III (AT-III) activity/assay, 97–98, 98*t*
- Antithyroglobulin antibody, 99, 99*t*
- Antithyroid drugs  
thyroid-stimulating hormone and, 887, 892  
total thyroxine (T<sub>4</sub>) and, 895  
WBC levels and, 976
- Antithyroid microsomal antibody, 100, 100*t*
- Antithyroid peroxidase antibody (anti-TPO), 100, 100*t*
- Antivinculin antibodies, 316
- Aortic artery pressure, 202–203*t*
- APC (activated protein C) resistance test, 397
- Apnea, sleep studies and, 827
- Apo A-I, 101–103
- Apo A-I/Apo B ratio, 101
- Apolipoprotein A (apo A), 102
- Apolipoprotein B (Apo B), 14, 101–102
- Apolipoprotein E (apo E), 102
- Apolipoproteins, 101–104, 103–104*t*
- Apt test, 105
- Argentaffin-staining (enteroendocrine) cells, 528
- Arginine stimulation, 460
- Arginine test, 477–478
- ARMD risk analysis, 20
- Arrhythmias, malignant ventricular, signal-averaged ECG and, 344
- Arsenic, 620, 620*t*
- Arsenicals, WBC levels and, 976
- Arterial blood gases (ABGs), 106–113, 108*t*, 112–113*t*
- Arterial blood oxygen saturation (SaO<sub>2</sub>), 658
- Arterial plethysmography, 713–714
- Arterial vascular balloon dilation, 114
- Arteriography (angiography), 114–117, 713  
bronchial, 757  
coronary, 201–206, 202–203*t*, 204*f*  
fluorescein, 424–425, 425*t*  
kidney, 117  
of peripheral vascular system, 116  
pulmonary, 757–758
- Arteriosclerotic peripheral vascular occlusion, 713
- Arthritis  
inflammatory, 118  
rheumatoid, 70, 790  
tests for, 1011
- Arthrocentesis, with synovial fluid analysis, 118–121, 121*t*
- Arthroscopy, 122–124, 122*f*, 124*t*
- Ascaris* spp. (hookworm), in stool, 847
- Ascitic fluid cytology, 671–674, 672*f*, 673–674*t*
- Ascorbic acid, triglycerides and, 908
- Asparaginase, triglycerides and, 908
- Aspartate aminotransferase (AST), 125–127, 126–127*t*
- Aspartate transaminase, 350
- Aspermia, 801
- Aspiration, 118  
bone marrow, 160–164, 160*t*, 163*f*, 164*t*  
pericardial, 689  
scan, 441–442  
suprapubic, of urine, 948
- Aspirin, 709  
in Bence Jones protein, 134  
free thyroxine (T<sub>4</sub>) and, 895  
thromboelastography and, 874–875  
thyroid-stimulating hormone and, 887  
WBC levels and, 976

## 1024 index

- Aspirin resistance tests, 709–711, 711*t*  
Atherectomy, of coronary arterial plaques, 203  
AT-III (antithrombin III) activity/assay, 97–98, 98*t*  
Atrial natriuretic peptide (ANP), 636–637, 637*t*  
Audiovisual recordings, 826  
Auditory brainstem-evoked potentials (ABEPs), 394  
    prolonged latency for, 395  
Australian antigen, 494  
Autoantibody(ies)  
    insulin, 327, 327*t*  
    panel, diabetes mellitus, 327, 327*t*  
Autoimmune diseases  
    antinuclear antibody in, 46*t*, 86  
    liver, 81  
Autoimmune endocrine diseases,  
    21-hydroxylase antibodies in, 530  
Autoimmune hepatitis, 81  
Autoimmune thyroiditis, antithyroglobulin antibody in, 99  
Autoimmune-mediated disease, antiparietal cell antibody in, 90  
Autologous blood, for transfusions, 154  
Autopsy, virtual, 274  
Autotransfusions, 153  
Azotemia, total blood volume and, 900
- B**  
B vitamins, 426  
B<sub>12</sub> deficiency, IF blocking antibody in, 545  
Bacille Calmette-Guérin (BCG)  
    immunization, 916  
*Bacillus anthracis*, in anthrax, 141, 142*t*  
BACTEC method, for tuberculosis culture, 917  
Bacteremia, blood cultures in, 147  
Bacteria, peritoneal fluid and, 671–672  
Bacterial cultures, 119  
Bad cholesterol, 566  
Band cells, 975  
Barbiturates, 853*t*  
    testosterone and, 862  
    vitamin D and, 972  
    WBC levels and, 976  
Barium  
    in abdominal ultrasound, 2  
    in bone mineral density, 158  
    column, defects in, 131  
    reflux of, 128, 131  
    for small bowel follow-through, 830  
    for upper gastrointestinal x-ray study, 922  
Barium enema (BE), 128–130, 130*t*  
Barium swallow, 131–132, 132*t*  
Base excess, 106, 110  
Basic metabolic panel (BMP), 1011  
Basic metabolic profile (BMP), 133  
Basopenia, 977–978*t*  
Basophilia, 977–978*t*  
Basophils, 975  
B-cells, 975  
BCR-ABL fusion gene, 195–196*t*  
Bence Jones protein, 134–135, 135*t*  
Benign peptic ulceration, 922  
Benign prostatic hypertrophy (BPH), PSA and, 743  
Benzodiazepines, 853*t*  
Bernstein test, 378  
Beta<sub>2</sub>-microglobulin, 195–196*t*, 628–629, 629*t*  
Beta-blockers, thromboelastography and, 874–875  
Beta-D-glucan, 428–429, 429*t*  
11 beta-prostaglandin F(2) alpha, 136, 136*t*  
Bethesda System, for Pap smear, 667  
Bethesda System for Reporting Thyroid Cytopathology, 881, 881*t*  
Bicarbonate (HCO<sub>3</sub>), 106–109, 108*t*, 197  
Bilateral adrenal nodular hyperplasia, aldosterone in, 25  
Bile, 137  
Bile ducts, 364, 366*f*  
    pancreatic enzymes and, 662  
    ultrasound of, 1–2, 4  
Biliary cirrhosis, primary  
    antimitochondrial antibody in, 82  
Bilirubin, 137–140, 139–140*t*  
    metabolism and excretion of, 137, 138*f*  
    urinary, 935–936, 939  
Biophysical profile, 406–408  
    fetal, 406–408  
Biopsy  
    bone marrow, 136, 160–164, 160*t*, 164*t*  
    cervical, 224–226  
    chorionic villus, 241–243, 242*f*  
    cone, 224–226  
    endocervical, 224  
    endometrial, 362–363  
    liquid, 569–570  
    liver, 571–573, 571*f*  
    lung, 584–587  
        closed technique of, 584  
        open method of, 584, 586  
        percutaneous, 843  
        percutaneous needle biopsy as, 586  
        thoroscopic, 586

- Biopsy (*Continued*)
- transbronchial, 585, 585*f*
  - transbronchial brushing, 586
  - transbronchial needle aspiration, 585, 585*f*
  - pleural, 715–716
  - punch, 224
  - renal, 776–778, 777*f*
  - sentinel lymph node, 805–806
  - simple cervical, 224
  - skin, 824
  - stereotactic, of breast, 606, 607*f*
- Bioterrorism infectious agents testing, 141–144, 142*t*
- Biotinidase deficiency, newborn metabolic screening and, 643
- Bismuth (Pepto-Bismol), urea breath test and, 925
- Bitewing view, 322
- Bladder
- cystography, 307
  - cystometry, 950
- Bladder cancer, genomics of, 456–457
- Bladder cancer markers, 145–146
- Bladder tumor antigen (BTA), 145–146
- Bleeding disorders, fibrinogen for, 422
- Bleeding time (BT), 709–711, 711*t*
- Blood
- bone turnover markers in, 171
  - in cerebrospinal fluid, 577
  - chloride in, 233–234, 234*t*
  - cortisol in, 290–292, 292*t*
  - creatinine in, 303
  - D-xylose in, 983
  - leucine aminopeptidase in, 561, 561*t*
  - mean platelet volume for, 712
  - metrapone test in, 17
  - neuron-specific enolase in, 638
  - osmolality of, 654–655, 655*t*
  - parathyroid hormone for, 675
  - partial thromboplastin time for, 681
  - parvovirus B19 antibody for, 684
  - pepsinogen for, 688
  - per cubic milliliter of, 706
  - pheochromocytoma suppression and provocative testing for, 692
  - phosphate for, 694
  - phospholipase A2 receptor antibodies for, 696
  - PI-linked antigen for, 697
  - placental growth factor for, 698
  - plasminogen activator inhibitor 1 for, 700
  - plasminogen for, 699
  - platelet aggregation test for, 702
  - platelet antibody detection for, 704
  - platelet count for, 706
  - platelet function assay and, 709
  - potassium for, 724–726, 726*t*
- Blood (*Continued*)
- prealbumin for, 728
  - pregnancy-associated plasma protein-A for, 730
  - progesterone assay for, 733
  - prolactin levels in, 737
  - prostate-specific antigen and, 743–745
  - protein C and S for, 751
  - prothrombin time for, 753
  - rabies-neutralizing antibody test in, 769
  - red blood cell count in, 770–771, 771*t*
  - red blood cell indices for, 772–775
  - renin assay in, 783–787
  - reticulocyte count in, 788
  - rheumatoid factor in, 790
  - ribosome P antibodies in, 792
  - rubella antibody test in, 793
  - SARS viral testing of, 797
  - serotonin in, 807
  - sexual assault testing of, 809
  - sickle cell screen of, 820
  - sodium in, 835–836, 836*t*
  - Streptococcus* serologic testing of, 850
  - substance abuse testing of, 852
  - syphilis detection test of, 859
  - testosterone in, 861
  - uric acid in, 927–929
  - Zika virus in, 985
- Blood alcohol, 391–392
- Blood antibody screening, 289
- Blood chromosome analysis, 244–245, 245*t*
- Blood clotting, platelet activity and, 706
- Blood crossmatching, blood typing and, 154
- Blood culture and sensitivity, 147–148
- Blood donors, anti-HTLV-I/II antibodies in, 525
- Blood EtOH, 391–392
- for alcohol levels, 391–392
- Blood folate, 426
- Blood gases, 106–113, 108*t*, 112–113*t*
- Blood group microarray testing, 152–154
- Blood indices, 772–775, 773*t*
- Blood sampling, fetal, fetoscopy for, 420*f*
- Blood smear, 149–151, 150*b*
- Blood spot, 642
- Blood sugar, 462–464, 464*t*
- Blood tests/testing, 1000–1005
- allergy, 31–32, 32*t*
  - for amino acids, 45
  - in *Helicobacter pylori* testing, 482
- Blood transfusions, 639
- blood typing and, 153
  - Coombs test and, 287
  - neutrophil antibody screen, 639
  - serum iron in, 547
- Blood typing, 152–154, 152*t*
- Blood urea nitrogen (BUN), 155–156, 156*t*

## 1026 index

Blood-clotting factors, 248–250, 248–250*t*, 252*t*

Blood-clotting mechanism, fibrinogen and, 422

Body movements, fetal, 407

Body plethysmography, 761–762

Bone

- alkaline phosphatase in, 29
- aluminum in, 43
- PET scan and, 721
- tests for, 1011

Bone collagen equivalents (BCEs), 168–171, 171*t*

Bone densitometry, 157–159

Bone Gla protein (BGP), 170

Bone marrow

- biopsy, 160–164, 160*t*, 164*t*
- in SMCD, 136
- erythropoietin and, 375
- transplantation, veno-occlusive disease after, 700

Bone mineral content (BMC), 157–159

Bone mineral density (BMD), 157–159

- testing of, 158

Bone scan, 165–167, 167*t*

Bone turnover markers (BTMs), 168–171, 171*t*

Bone x-ray, 172–173, 173*t*

Bone-specific alkaline phosphatase (BSAP), 170

- serum, 169

*Borrelia burgdorferi*, 595

Botulism infection, bioterrorism infectious agents testing for, 141, 142*t*

Bowel perforation, amylase in, 56

*BRAF*, 588

- mutation analysis, 256–257
- testing for, 588–589

*BRAF* V600 mutations, 195–196*t*

Brain

- aluminum in, 43
- computed tomography of, 278–280, 280*t*
- electroencephalography of, 347
- evoked potentials of, 394
- MRI of, 599
- oxygen testing and monitoring, 658
- PET for, 719

Brain natriuretic peptide (BNP), 636–637, 637*t*

Brain scan, 174–175, 175*t*

Brain trauma, 719–720

BRCA (breast cancer) genes, 445–452, 446*t*

Breast cancer

- ductoscopy in, 336, 337*f*
- estrogen receptor assay in, 389
- genetic testing in, 445–447, 446*t*, 454

Breast cancer (*Continued*)

- mammography for, 606
- detection rate for, 606
- progesterone receptor assay for, 735
- sentinel lymph node biopsy for, 805
- tumor analysis, 176–178

Breast ductal lavage, 179–180, 336

Breast MRI, 600

Breast sonogram, 181–182, 182*t*

Breast tomography, 606–607

Breast ultrasonography, 181–182, 182*t*

Breath, for drug testing, 852

Breath test, 391

- for *Helicobacter pylori*, 483

Breathing movements, fetal, 406–407

Bronchial angiography, 757

Bronchial provocation studies, 763

Bronchoalveolar lavage, for SARS viral testing, 797

Bronchoscopy, 183–186, 185*f*, 186*t*, 843

- virtual, 281

*Brucella* sp., 142*t*, 143

Brucellosis, bioterrorism infectious agents testing for, 142*t*, 143

Brushing, transbronchial, 586

BT (bleeding time), 709–711, 711*t*

BUN/creatinine ratio, 155

Burkitt lymphoma

- EA-R and, 368
- EBV and, 368

Burr cells, in red blood cells, microscopic examination of, 150*b*

## C

C1 complement, deficiency, 267*t*

C1-INH complement, deficiency, 267*t*

C1q complement, deficiency, 267*t*

C1r complement, deficiency, 267*t*

C1s complement, deficiency, 267*t*

C2 complement, deficiency, 267*t*

C3 complement, 266–268, 267–268*t*

C4 complement, 266–268, 267–268*t*

C5 complement, deficiency, 267*t*

C6 complement, deficiency, 267*t*

C7 complement, deficiency, 267*t*

C8 complement, deficiency, 267*t*

C9 complement, deficiency, 267*t*

CA15-3/CA27.29, 195–196*t*

CA19-9, 195–196*t*

CA-125, 195–196*t*

Caffeine

- catecholamines and, 956
- vanillylmandelic acid and, 955

Calcifications, in bone x-ray, 172

Calcitonin, 187–188, 188*t*, 195–196*t*

- Calcium (Ca), 189–191, 189*t*, 191*t*, 252*t*  
 parathyroid hormone level and, 675  
 with potassium, and magnesium, 597  
 vitamin D and, 971
- Calcium infusion test, 187, 440
- Calculus analysis, renal, 945–946
- Caloric study, 192–193
- Calprotectin, fecal, 400–401, 401*t*
- Campylobacter*-like organism (CLO) test, 482–484
- Cancer, 677
- bladder
    - genomics of, 456–457
    - markers, 145–146
    - superficial, recurrence rate for, 145
  - breast, 445–452, 446*t*
    - ductoscopy in, 336, 337*f*
    - estrogen receptor assay in, 389
    - genetic testing in, 445–447, 454
    - mammography in, 606
    - progesterone receptor assay for, 735
    - sentinel lymph node biopsy for, 805
    - tumor analysis, 176–178
  - cervical, Pap test for, 520, 667
  - colon, 445–452, 446*t*
    - genetic testing of, 447, 455
    - tumor analysis, 256–257
  - cytologic test for, 667–670, 670*t*
  - lung
    - genomic testing of, 455
    - molecular testing for, 588–589
    - neuron-specific enolase in, 638
  - lymphoma and hematologic, genomics of, 456
  - ovarian, 445–452, 446*t*
    - genetic testing in, 445–447
  - prostate, 739
    - acid phosphatase in, 7
    - genomic testing of, 454–455
    - metastatic, 739
    - PSA in, 743–745, 745*t*
  - screening, stool, 844–846
  - studies for, 988
  - thyroid, 445–452, 446*t*, 879
    - genetic testing of, 449–450, 455–456
  - tumor markers, 194–196, 195–196*t*
- Candida*, 813*t*
- Capsule endoscopy, 382
- Captopril
  - renal scan, 779–782
  - test, 785
- Carbon (<sup>13</sup>C)-labeled urea, 925
- Carbon dioxide (CO<sub>2</sub>) combining power, 197–198, 198*t*
- Carbon dioxide (CO<sub>2</sub>) content, 197–198, 198*t*
- Carbon dioxide (CO<sub>2</sub>) monitor, 826
- Carbon monoxide, 199–200
  - poisoning, 199
- Carboxyhemoglobin (COHb), 199–200
- Carcinoembryonic antigen (CEA), 145–146, 195–196*t*
  - peritoneal fluid and, 671–672
  - in pleural fluid analysis, 867
- Carcinoid nuclear scan, 652–653
- Carcinoid syndrome, 807
- Carcinoid tumors, 807
  - 5-HIAA in, 528
- Cardiac arrhythmias, and magnesium, 597
- Cardiac catheterization, 201–206, 202–203*t*, 204*f*
- Cardiac conduction slowing, 597
- Cardiac echo, 339–341
- Cardiac enzyme, creatinine kinase as, 299
- Cardiac flow studies, 207–209
- Cardiac genetic testing, 450
- Cardiac index (CI), 202–203*t*
- Cardiac injury, tests for, 1011
- Cardiac mapping, 359–361
- Cardiac monitoring, pressures and volumes used in, 202–203*t*
- Cardiac nuclear scan, 207–209
- Cardiac output (CO), 202–203*t*
- Cardiac perfusion scan, 207
- Cardiac risk factors
  - lipid profile, 565–568
  - lipoprotein-associated phospholipase A<sub>2</sub>, 564
- Cardiac scan, 207–209
- Cardiac stress testing, 210–213, 339–340
- Cardiac tamponade, 689
- Cardiac tumors, ultrasound for, 339
- Cardiac-specific troponin I (cTnI), 912–914. *See also* Troponins
- Cardiac-specific troponin T (cTnT), 912–914. *See also* Troponins
- Cardiology, PET scan and, 720
- Cardiomyopathy
  - antimyocardial antibody in, 83
  - hypertrophic, 450
- Cardiovascular disease, 445–452, 446*t*
  - fibrinogen and, 422
  - lipoprotein-associated phospholipase A<sub>2</sub> in, 564
- Cardiovascular system, tests for, 989
- Carotid artery duplex scanning, 214–215
- Carotid intima-media thickness (CIMT), 214
- Carotid ultrasound, 214–215
- Carrier couple, 447
- Cascara sagrada, urine color and, 932*t*
- Casts, urinary, 936–937
  - cellular, 936
  - epithelial, 943

## 1028 index

- Casts, urinary (*Continued*)  
fatty, 936, 943  
granular, 936, 943  
hyaline, 936, 943  
red blood cell, 937  
tubular (epithelial), 937  
waxy, 936, 943  
white blood cell, 937
- Catecholaminergic polymorphic ventricular tachycardia (CPVT), 450
- Catecholamines, 954–957
- Cathepsin D, 176–178
- Catheterization  
cardiac, 201–206, 202–203*t*, 204*f*  
ureteral, by cystoscopy, 310*f*
- C-cell hyperplasia, 187
- CCP IgG anti-CCP, 70, 70*t*
- CD4 count, antiretroviral therapy based on, 504*t*
- CD4 marker, 219–221, 219*t*, 221*t*
- CD4 percentage, 219–221, 219*t*, 221*t*
- CD4/CD8 ratio, 219–221, 219*t*, 221*t*
- CD20, as cancer tumor marker, 195–196*t*
- Celiac disease, 458
- Cell counts  
in pleural fluid analysis, 865  
in synovial fluid, 119
- Cell culture drug resistance testing (CCDRT), 216
- Cell size, 772
- Cell surface immunophenotyping, 219–221, 219*t*, 221*t*
- Cell-free DNA in maternal blood, 217–218
- Cell-free maternal DNA testing, 217–218
- Cellular casts, 936
- Cellularity, estimation of, in bone marrow biopsy, 161–162
- Central sleep apnea, 827
- Central venous pressure, 202–203*t*
- Centromere antibody, 67
- Cephalosporins, cold agglutinins and, 253
- Ceramides, 222–223, 222*t*
- Cerebral arteriography, brain computed tomography and, 278
- Cerebral blood flow, 174–175, 175*t*  
electrophysiologic effects of, EEG and, 347
- Cerebral lesions, EEG and, 347
- Cerebrospinal fluid (CSF)  
analysis of  
contraindications to, 580  
potential complications of, 580  
prealbumin, 728  
blood in, 577  
cells in, 578  
chloride in, 578  
color of, 577  
C-reactive protein in, 579
- Cerebrospinal fluid (CSF) (*Continued*)  
culture and sensitivity of, 578  
cytology of, 579  
examination, 576–583  
glucose in, 578  
glutamine, 579  
lactic acid, 579  
lactic dehydrogenase in, 579  
lumbar puncture and, 576–583, 581*f*  
pressure of, 577  
protein in, 578  
serology for syphilis in, 579  
*Streptococcus* serologic testing of, 850  
tumor markers for, 579  
Zika virus in, 985
- Cervical biopsy, 224–226
- Cervical block, 225
- Cervical cancer, Pap test for, 520, 667
- Cervical culture, 815  
for *Chlamydia*, 231  
for herpes simplex, 499
- Cervical mucus sperm penetration test, 822–823
- Cervical screening, recommendations for, 521*t*
- Cervical spinal x-ray studies, 839
- Cervical spine, MRI of, 601
- Cervicography, 263–264
- Cervix  
cytology of, 667–668  
Pap test for, 667–670, 669*f*, 670*t*  
cf DNA screening, 217–218
- CH 50, complement, 266
- Chem 7, 133
- Chem 12, 270–271, 270*t*
- Chemical stress testing, 210
- Chemistry panel, 133, 270–271, 270*t*
- Chemistry screen, 133
- Chemosensitivity assay, 216
- Chemotherapeutic agents, WBC levels and, 976
- Chemotherapy  
cell culture drug resistance testing, 216  
high-dose, 700
- Chest  
computed tomography of, 281–283  
x-ray of, 227–229, 229*t*
- Chest impedance, 826
- Chest pain, ischemia-modified albumin in, 550
- Chest x-ray, lung scan and, 590
- CHF peptides, 636–637, 637*t*
- Chlamydia*, 230–232, 813*t*  
in sexually transmitted disease testing, 813–814
- Chlamydia trachomatis*, 813*t*
- Chlamydia pneumoniae*, 230



- Chlamydomytila psittaci*, 230  
*Chlamydomytila trachomatis*, 230  
 Chloral hydrate, catecholamines and, 956  
 Chloride (Cl)  
   blood, 233–234, 234*t*  
   in cerebrospinal fluid, 578  
   values in children, 85*t*  
 Chloroform, WBC levels and, 976  
 Chloroquine (Aralen), urine color and, 932*t*  
 Chlorpromazine, in metyrapone test, 17  
 Cholescintigraphy, 431–432  
 Cholestatic drugs, urobilinogen and, 939  
 Cholesterol, 235–237, 236–237*t*  
 Cholesterol-to-HDL ratio, 235, 236*t*  
 Cholestyramine  
   triglycerides and, 908  
   vitamin D and, 972  
 Cholinesterase, 238–240, 240*t*  
   RBC, 238–240, 240*t*  
 Chorion frondosum, 241  
 Chorionic villus biopsy (CVB), 241–243, 242*f*  
 Chorionic villus sampling (CVS), 241–243, 242*f*  
 Christmas factor, 249–250*t*, 252*t*  
 Chromium, aluminum and, 43  
 Chromogranin A, 195–196*t*, 807–808, 808*t*  
 Chromophobic adenoma, prolactin and, 737  
 Chromosomal aberrations, 51  
 Chromosomal abnormalities, women at  
   high risk for having baby with, 217, 218*b*  
 Chromosome karyotype, 244–245, 245*t*  
 Chromosome studies, 244–245, 245*t*  
 Chromosomes 3, 17, and 9p21, 195–196*t*  
 Chronic active hepatitis (CAH),  
   autoimmune, antismooth muscle  
   antibodies in, 92  
 Chronic disease, anemia of, 906  
 Chronic GM2, 447  
 Chronic hepatitis, anti-LKM-1 antibodies  
   in, 81  
 Chronic illness, iron level in, 548  
 Chronic inflammatory pattern, protein  
   electrophoresis and, 748*t*  
 Chronic iron overload, 547  
 Chronic kidney disease (CKD), 302  
   neutrophil gelatinase-associated lipocalin  
   for, 640  
 Chronic primary adrenal insufficiency,  
   21-hydroxylase antibodies in, 530  
 Chronic renal diseases, water load test in, 72  
 Chymotrypsin, 662–664, 664*t*  
 Circadian rhythm, in amino acids, 45  
 Cisternogram, 174–175, 175*t*  
 CK-BB, 298  
 C-kit/CD117, 195–196*t*  
 CK-MB, 297  
 CK-MM, 298–299  
 Cl (chloride)  
   blood, 233–234, 234*t*  
   in cerebrospinal fluid, 578  
 Claus breast cancer risk models, 179  
 Clean-catch urine collection, 939, 948  
 CLO test, for *Helicobacter pylori*, 482  
 Clofibrate  
   thromboclastography and, 874–875  
   total thyroxine ( $T_4$ ) and, 895  
   triglycerides and, 908  
 Clomiphene (Clomid) Challenge Test  
   (CCCT), 593–594  
 Clonidine, 692  
   catecholamines and, 956  
   vanillylmandelic acid and, 955  
 Clonidine suppression test (CST), 692–693  
 Clostridial toxin assay, 246–247, 247*t*  
*Clostridium botulinum*, in botulism, 141, 142*t*  
*Clostridium difficile*, 246–247, 247*t*  
*Clostridium difficile*-associated disease  
   (CDAD), 246  
 CMV-specific IgG avidity, 317  
 CO (cardiac output), 202–203*t*  
 Coagulating factors, 248–250, 248–250*t*, 252*t*  
 Coagulating factors concentration,  
   248–250, 248–250*t*, 252*t*  
 Coagulation, antithrombin III in, 97  
 Coagulation factor, PTT test, 681  
 Coagulation factor inhibitors, 250  
 Coagulation panel, 255, 255*t*  
 Coagulation profile, 255, 255*t*  
 Coagulation screening, 1011  
 Cocaine, 853*t*  
 Coccygeal x-ray studies, 839–840  
 Codocytes, in red blood cells, microscopic  
   examination of, 150*b*  
 COHb, 199–200  
 Coitus interruptus, 804  
 Cold agglutinin disease, 253  
 Cold agglutinins, 253–254, 254*t*  
 Cold water, caloric study and, 192  
 Colestipol, triglycerides and, 908  
 Colon, barium enema in, 128  
 Colon cancer, 445–452, 446*t*  
   genetic testing of, 447, 455  
   tumor analysis, 256–257  
 Colon Recurrence Score test, 455  
 Colonoscopy, 258–261  
   virtual, 273  
 Color Doppler echocardiography, heart, 339  
 Color Doppler ultrasound (CDU), 214, 960  
   scrotal, 799

## 1030 index

- Color flow imaging, useful application of, [339](#)
- Colorectal cancer
  - fecal occult blood test for, [844](#)
  - hereditary nonpolyposis, [256](#)
- ColoVantage, [625](#), [625t](#)
- Colpography, video, [263–264](#)
- Colposcope, [810](#)
- Colposcopy, [262–265](#), [262f](#)
  - video, [263–264](#)
- Coma, tests for, [1011](#)
- Committee on Thrombolytic Agents (CTA), [699](#)
- Common associative diseases with thrombocytosis, [706](#)
- Complement assay, [266–268](#), [267–268t](#)
- Complement levels, synovial fluid in, [119](#)
- Complete blood count (CBC), [269](#), [772](#)
- Complete deficiency, of PAI-1, [700](#)
- Comprehensive metabolic panel (CMP), [270–271](#), [270t](#), [1012](#)
- Comprehensive pulmonary function study, [760](#)
- Computed tomography (CT)
  - of abdomen and pelvis, [272–277](#), [277t](#)
  - of brain, [278–280](#), [280t](#)
  - of chest, [281–283](#)
  - four-dimensional, [281–282](#)
  - of heart, [284–286](#)
  - of parathyroid glands, [281–282](#)
  - of pelvis, [272–277](#), [277t](#)
  - PET scan and, [719](#)
  - quantitative, in bone mineral density, [157](#)
  - single-photon emission, [433](#)
    - cardiac, [207](#)
    - of liver, [574](#), [581f](#)
- Computer-assisted semen analysis, [802](#)
- Conduction velocity, [355](#)
- Cone biopsy, [224–226](#)
- Congenital adrenal hyperplasia (CAH)
  - adrenal steroid precursors in, [18](#)
  - newborn metabolic screening and, [643](#)
- Congenital blood disorders, fetal blood samples for, [419](#)
- Congenital deficiencies, of vitamin K-dependent proteins, [751](#)
- Congenital hypothyroidism, newborn metabolic screening and, [643](#)
- Congenital immune deficiencies, immunoglobulins in, [537](#)
- Congenital syndromes, [702](#)
- Congestive heart failure
  - natriuretic peptide, [636](#)
  - total blood volume for, [900](#)
- Conization, [224](#)
- Conjugated (direct) bilirubin, [137](#)
- Conn syndrome, aldosterone in, [25](#)
- Connecting peptide insulin, [293–294](#), [294t](#)
- Contraction stress test (CST), [409–411](#)
- Contrast media, catecholamines and, [956](#)
- Conventional nuclear bone scans, [721](#)
- Conventional Pap smear testing (CPT), [667–668](#)
- Coombs test
  - direct, [287–288](#)
  - indirect, [289](#)
- Core window, [494](#)
- Coronary angiography, [201–206](#), [202–203t](#), [204f](#)
- Coronary angioplasty, percutaneous transluminal, [203](#)
- Coronary arteries, intravascular ultrasound in, [543](#)
- Coronary artery disease, high-density lipoproteins in, [565](#)
- Coronary atherectomy, intravascular ultrasound in, [543](#)
- Coronary atherosclerosis, intravascular ultrasound in, [543](#)
- Coronary calcium score, [284–286](#), [285t](#)
- Coronary CT angiography, [284–286](#)
- Coronary heart disease (CHD), cholesterol-to-HDL ratio for, [235](#), [236t](#)
- Coronary or peripheral artery disease (CPAD), apo A in, [102](#)
- Coronary stent placement, intravascular ultrasound in, [543](#)
- Coronavirus (CoV), [797](#)
- Corticosteroids, vitamin D and, [972](#)
- Corticotropin, [11–12](#), [12t](#)
- Corticotropin-releasing hormone (CRH), [11](#)
- Cortisol, [290–292](#), [292t](#)
  - in 17-OHCS, [526](#)
  - in adrenocorticotropic hormone, [11](#)
- Cortisol stimulation test, [13–15](#)
- Cortisol suppression test, [324–326](#), [326t](#)
- Cortrosyn, [13](#)
- Cosyntropin, adrenocorticotropic hormone stimulation test with, [13–15](#)
- Cotinine, [646–648](#), [646t](#)
- Coumadin, thromboelastography and, [874–875](#)
- COX1 (cyclo-oxygenase-1), [710](#)
- C-peptide, [293–294](#), [294t](#)
- CPK (creatin phosphokinase), [297–300](#), [298f](#), [300t](#)
  - aldolase and, [23](#)
  - isoenzyme, [635](#)
- Cranial nerve, eighth, [192](#)
- C-reactive protein
  - in alpha defensin test, [39](#)
  - in cerebrospinal fluid, [579](#)
  - in lipoprotein-associated phospholipase A<sub>2</sub>, [564](#)

- C-reactive protein test, 295–296, 296*t*, 373
- Creatine kinase (CK), 297–300, 298*f*, 300*t*
- Creatine kinase-MB (CK-MB), troponins and, 912
- Creatine phosphokinase (CPK), 297–300, 298*f*, 300*t*
- aldolase and, 23
- isoenzyme, 635
- Creatinine, blood, 301–302, 302*t*
- Creatinine clearance (CC), 303–305, 305*t*
- CREST syndrome, anticentromere antibody test in, 67
- Crohn disease, lactoferrin, 555
- Crohn disease prognostic panel, 80, 80*t*
- Crossmatching, of blood, 154
- Cross-table lateral view, of abdomen, 650
- Cryoglobulin, 306
- Cryoglobulinemia, type I (monoclonal), 306
- Crystals, urinary, 936, 939, 943
- CT (computed tomography) arteriography, 273–274, 278
- CT (computed tomography) colonography, 273
- CT (computed tomography) nephrotomography, 274
- CT (computed tomography) urography, 274
- C-telopeptide (CTx), 168–169
- C-type natriuretic peptide (CNP), 636–637, 637*t*
- Culture, 814
- anal canal, 815
- bacterial, 119
- cervical, 815
- for *Chlamydia*, 231
- for herpes simplex, 499
- fungal, 119, 428–429, 429*t*
- Gram stain and bacteriologic, 866
- for herpes simplex, 498
- cervical, 499
- urethral, 499
- for *Mycobacterium tuberculosis* and fungus, 866
- nose, 872–873
- oropharyngeal, 815
- rectal, 816*f*
- sexually transmitted disease, 814
- sputum, 841–842
- stool, 847–848
- throat, 872–873
- tuberculosis, 917–918
- urethral, 814, 816, 817*f*
- urine, 816, 947–948
- of vaginal secretions, 814
- viral, 965–966, 965*t*
- wound, 981–982
- Culture and sensitivity (C&S)
- of cerebrospinal fluid, 578
- sputum, 841–842
- stool, 847–848
- urine, 947–948
- wound, 981–982
- Curettage, endocervical, 224
- Cushing syndrome
- 17-OHCS in, 526
- adrenocorticotrophic hormone in, 11
- adrenocorticotrophic hormone stimulation test with cosyntropin in, 13, 15
- cortisol and, 290
- cortisol suppression test for, 324
- metyrapone test in, 16
- Cutaneous anthrax, 141
- Cutaneous immunofluorescence
- biopsy, 824
- Cyanocobalamin, 969–970
- Cyclic citrullinated peptide antibody, 70, 70*t*
- Cyclo-oxygenase-1 (COX1), 710
- CYP450 (cytochrome P450) system, 334
- Cystatin C, 302
- Cystic fibrosis (CF), 445–452, 446*t*, 662, 857
- genetic testing of, 448
- newborn metabolic screening and, 644
- Cystic fibrosis transmembrane conductance regulator (CFTR) gene, 448
- Cystography, 307–308
- Cystometrogram (CMG), 949–952
- Cystometry, 949–952
- Cystoscopy, 309–313, 309–310*f*
- Cystourethrography, 307–308
- Cytochrome P450 (CYP450) system, 334
- Cytogenetics, 244–245, 245*t*
- Cytokeratin fragment 21-1, 195–196*t*
- Cytokines, 314–315
- Cytolethal distending toxin B (CdtB), 316
- Cytology
- in bladder cancer, 145
- in pleural fluid analysis, 867
- Cytomegalovirus (CMV), 317–318
- TORCH test for, 899
- Cytoplasmic ANCA (c-ANCA), 84
- ## D
- Danazol
- free thyroxine (T<sub>4</sub>) and, 895
- thyroxine-binding globulin and, 897
- Dane particle, 494
- DaT scan, 174–175, 175*t*, 679
- D-dimer test, 319–320, 320*t*, 877
- De Ritis ratio, 21
- Dead space, 761
- Decubitus images, in chest x-ray, 227

## 1032 index

- Deep vein thrombosis (DVT)
  - D-dimer test, 319
  - vascular ultrasound studies for, 960
- Dehydroepiandrosterone (DHEA), 18–19,  
18–19*t*, 861
- Dehydroepiandrosterone sulfate (DHEA S),  
18–19, 18–19*t*
- 11-Dehydro-thromboxane B2 (11-dTXB2),  
709–711, 711*t*
- Delta gap, 654–655
- Delta hepatitis, 497
- Delta-aminolevulinic acid (d-ALA), 321,  
321*t*
- Dental radiography, 322–323
- Dental x-ray, 322–323
- 11-Deoxycortisol, 18–19, 18–19*t*
- Dermatitis, allergy blood testing in, 31
- Dermatographism, allergy blood testing  
in, 31
- Dermographism, allergy skin testing in,  
33–34
- Des-gamma-carboxy prothrombin,  
195–196*t*
- Dexamethasone, testosterone and, 862
- Dexamethasone suppression test, 291,  
324–326
- Dextran, thromboelastography and,  
874–875
- Diabetes insipidus (DI)
  - antidiuretic hormone in, 71
  - nephrogenic, 934–935
- Diabetes mellitus (DM)
  - albumin and, 626
  - arginine stimulation and, 461*t*
  - autoantibody panel, 327, 327*t*
  - fasting blood sugar in, 462–464
  - gestational, 465
  - glucagon in, 462–464
  - glucose monitoring in, 462–464
  - glucose tolerance test in, 467–470, 468*f*
  - glycosylated hemoglobin in, 471–473,  
472–473*t*
  - human placental lactogen in, 523
  - management of, tests for, 1012
  - microalbumin and, 626
  - postprandial glucose, 465–466
- Diabetic control index, 471–473, 472–473*t*
- Diabetic nephropathy, microalbumin and,  
626
- Diagnostic arthroscopy, 122–123
- Diamox (acetazolamide), 174
- Diarrhea, D-xylose absorption test for, 983
- Diarrhea-predominant irritable bowel  
syndrome (D-IBS), 316
- Diathesis, bleeding, 709
- Diatrizoate, in barium enema, 128
- Dibucaine, 238
- Dibucaine inhibition number, 238
- Dibucaine inhibition test, for cholinesterase,  
238
- Dicumarol, coumarin derivatives, 754
- Dietary deficiency, of folate, in  
homocysteine, 513
- Dietary phosphorus, 694
- Diethylstilbestrol, testosterone and, 862
- Differential count, 269
  - white blood cell count and, 974–978,  
974*t*
- Differentiation, of total bilirubin, 138–139
- Digital mammography, 606
- Digital subtraction angiography (DSA), 114
- Digoxin, 333*t*
  - testosterone and, 862
- Dihydrotestosterone (DHT), 861, 861*t*,  
863*t*
- 1,25-Dihydroxyvitamin D (1,25(OH)<sub>2</sub>D),  
971–973
- Dilated cardiomyopathy (DCM), 450
- 2,3-diphosphoglycerate, 328–329, 329*t*
- Dipyridamole (Persantine), 211
- Direct antiglobulin test (DAT), 287–288
- Direct bilirubin, 137
- Direct Coombs test, 287–288
- Direct immunofluorescence antibody test,  
824
- Direct lateral view, 610
- Direct viral testing, for HIV infection, 503
- Direct visualization, of fetus, 419, 420*f*
- Directed blood, for donation, 154
- Diseases
  - associated with complement deficiencies,  
266–267, 267*t*
  - organ panels and, 1011–1014
- DISIDA (diisopropyl iminodiacetic acid  
analogue) scanning, 431–432
- Disseminated intravascular coagulation  
(DIC), 319, 331*f*, 331*t*, 877
  - screening for, 330, 331*t*
  - tests for, 1012
  - thrombosis indicators fibrin monomers  
for, 877
- Disseminated neonatal herpes virus  
infections, 498
- Distal latency, 355
- Disulfiram, catecholamines and, 956
- Diuretic renal scan, 781
- Diuretics, WBC levels and, 976
- Diversion, urinary, 948
- DNA antibody, 74, 74*t*
- DNA evidence collection, for sexual assault  
testing, 811*b*
- DNA home test, of stool, 846
- DNA mismatch repair (MMR) genetic  
testing, 256–257
- DNA ploidy status, 176–178
- DNA stool sample test, 844

DNA testing, 450  
 for HPV, 520–522  
 stool sample for, 844–846

Dobutamine, 211

Docusate calcium (Doxidan, Surfak), urine  
 color and, 932*t*

Dopamine  
 thyroid-stimulating hormone and, 887  
 urinary, 954–957

Doppler echocardiography, color, 339

Doppler ultrasound  
 of placenta, 408  
 scrotal, 799  
 venous, 960–962

Double test, 612

Double-stimulated test, 477

Down syndrome, cell-free maternal DNA  
 testing for, 217

Downey test, 105

Doxorubicin (Adriamycin), urine color  
 and, 932*t*

Dressler syndrome, antimyocardial antibody  
 in, 83

Drug monitoring, 332–335, 333*t*  
 genotyping for, 334, 335*t*

Drug resistance testing  
 cell culture, 216  
 in HIV, 502

Drug response assay, 216

Drug screening, 852–854, 853*t*

Drug sensitivity testing  
 genotype, 332–335  
 in HIV, 502

Drug-induced thrombocytopenia, 704

Drugs  
 bilirubin and, 139  
 calcitonin level and, 188  
 cardiac stress testing and, 212  
 cholesterol level and, 236  
 cholinesterase and, 239  
 Coombs test and  
 direct, 287–288  
 indirect, 289  
 cortisol and, 291  
 C-peptide and, 294  
 C-reactive protein test and, 296  
 creatine kinase level and, 299  
 creatinine and, 302  
 estimated glomerular filtration rate and,  
 305  
 serum calcium level and, 190–191  
 serum chloride level and, 233  
 serum CO<sub>2</sub> level and, 197–198  
 uric acid and, 928  
 urinary proteins and, 938  
 urinary specific gravity and, 938  
 urine color and, 932*t*  
 urine pH and, 938

DSMA renal scan, 779–782

DTPA renal scan, 779–782

Dual-energy x-ray absorptiometry (DEXA)  
 scan, 157–159

Ductal lavage, breast, ductoscopy, 336

Ductoscopy, 336–338, 337*f*, 338*t*

Duodenal obstruction, amylase in, 56

Duodenoscope, fiberoptic, 365, 366*f*

Duplex scan  
 carotid artery, 214–215  
 venous, 960–962

Dwarfism, 474

D-xylose, 983

D-xylose absorption test, 983–984,  
 983–984*t*

Dynamic CT scanning, 273–274

## E

E<sub>1</sub> (estrone), 386*t*, 387

E<sub>2</sub> (estradiol), 386–388, 386*t*, 388*t*

E<sub>3</sub> excretion (estriol), 386*t*, 387

Early antigen (EA), 368

Early prostate cancer antigen (EPCA), 743

EBV (Epstein-Barr virus) antibody titer,  
 368–370, 369*t*

EBV nuclear antigen (EBNA), 368

Echinocytes, in red blood cells, microscopic  
 examination of, 150*b*

Echo stress testing, 210–213

Echocardiography, 339–341  
 color Doppler, 339  
 M-mode, 339  
 three-dimensional, 339  
 transthoracic (TEE), 340, 903–905,  
 904*f*  
 transthoracic, 339–341  
 two-dimensional, 339

Eclampsia, plasminogen and, 699

ECoG (electrocorticography), 347

Ecstasy (MDMA), 853*t*

Ectopic pregnancy, human chorionic  
 gonadotropin in, 515–516

Edematous states  
 aldosterone in, 26  
 water load test in, 72

Edwards syndrome, 217

Ejection fraction (EF), 202–203*t*

EL1, 660–661, 661*t*

Elastase-1, 660–661, 661*t*

Electrocardiogram, 342–346, 344*f*, 724  
 in Holter monitoring, 510, 510*f*

Electrocardiography (ECG), 342–346,  
 344*f*, 826

Electrocorticography (ECoG), 347

Electrodes, for sleep studies, 827

Electrodiagnostic tests, 1005

Electroencephalogram, 347–349

## 1034 index

- Electroencephalography (EEG), 347–349, 826  
  electrodes, 347  
  for seizure disorder, 349
- Electrolyte, 724
- Electrolyte panel, 1012
- Electrolytes test, sweat, 857–858
- Electromyography (EMG), 347–349, 826  
  of pelvic floor sphincter, 353–354  
  for fecal incontinence, 353–354  
  for urinary incontinence, 353
- Electromyoneurography, 355
- Electroneurography (ENG), 355–356
- Electronystagmography, 357–358
- Electrooculography, 357–358, 826
- Electrophoresis  
  of CSF protein, 578  
  in immunoglobulin quantification, 536  
  lipoprotein, 565–568  
  protein, 746–750, 748*t*, 750*t*  
  for AAT deficiency, 37  
  of urine, 134
- Electrophysiologic study (EPS), 359–361
- Electroretinography, 393–394
- Elevated fecal calprotectin, 400
- Elevated NSE concentrations, 638
- Elliptocytes, in red blood cells, microscopic examination of, 150*b*
- Emphysema, alpha<sub>1</sub>-antitrypsin in, 37
- Empyema, in pleural fluid analysis, 864–865
- End-diastolic left ventricular pressure, 202–203*t*
- End-diastolic volume (EDV), 202–203*t*
- Endobronchial ultrasound, 183–184
- Endocervical biopsy, 224
- Endocervical curettage, 224
- Endocervical mucus sample, 822
- Endocrine system, tests for, 990–991
- Endogenous anticoagulant protein C, 397
- Endometrial ablation, hysteroscopy in, 534
- Endometrial biopsy, 362–363
- Endometrium  
  proliferative-type, 362  
  secretory-type, 362
- Endomysial antibodies (EMA), 458–459, 458*t*
- Endoscopic retrograde  
  cholangiopancreatography (ERCP), 364–367, 366*f*, 665  
  for jaundice, 364
- Endoscopy, 1005
- Endourethral urologic ultrasound, 1
- Endourology, 309–313, 309–310*f*
- End-systolic volume (ESV), 202–203*t*
- Enema  
  barium, 128–130, 130*t*  
  small bowel, 830–832
- Enoxaparin, 77
- Enterobius*, 848
- Enterochromaffin (EC) cells, 807
- Enteroendocrine cells, 528
- Enteroscopy, 382
- Enzyme(s)  
  cardiac, 299  
  drug metabolism by, 334  
  pancreatic, 662–664, 664*t*
- Enzyme immunoassay (EIA), 859  
  in antiglomerular basement membrane antibodies, 79
- Eosinopenia, 977–978*t*
- Eosinophilia, 977–978*t*
- Eosinophils, 975
- Epidermal growth factor (EGFR), 588  
  testing for, 588–589
- Epilepsy, PET scan and, 719–720
- Epinephrine  
  catecholamines and, 956  
  urinary, 954–957  
  vanillylmandelic acid and, 955  
  WBC levels and, 976
- Epithelial casts, urinary, 943
- Epithelial cells  
  abnormality, 224  
  of salivary gland, 795
- EPS (electrophysiologic study), 359–361
- Epstein-Barr virus (EBV), 630  
  infections, 368  
  testing, 368–370, 369*t*  
  mononucleosis, 368
- ER assay, 389–390
- Erect abdominal image, of x-ray study, 650
- Erythema chronicum migrans, 595
- Erythema infectiosum, in parvovirus B19 antibody, 684
- Erythrocyte count, 770–771, 771*t*
- Erythrocyte fragility, 371–372, 372*t*
- Erythrocyte sedimentation rate (ESR), 373–374, 374*t*
- Erythrocytes  
  2,3-DPG in, 328–329, 329*t*  
  peripheral blood smear in, 149  
  porphyrin fractionation of, 717
- Erythrocytosis  
  hematocrit in, 485  
  hemoglobin in, 488
- Erythromycin, catecholamines and, 956
- Erythropoietin (EPO), 375–376, 376*t*
- Esophageal function studies, 377–381, 380*f*
- Esophageal manometry, 377–381, 380*f*
- Esophageal motility studies, 377–381, 380*f*
- Esophageal pH monitoring, 377, 380*f*
- Esophageal pH probe, 826
- Esophagogastroduodenoscopy (EGD), 382–385
- Esophagoscopy, virtual, 281
- Esophagus, barium swallow in, 131

- Estimated glomerular filtration rate (eGFR),  
303–305, 304–305*t*
- Estradiol (E<sub>2</sub>), 386–388, 386*t*, 388*t*
- Estradiol receptor, 389–390
- Estriol, human placental lactogen and, 523
- Estriol (E<sub>3</sub>) excretion, 386*t*, 387
- Estrogen fractions, 386–388, 386*t*, 388*t*
- Estrogen receptor assay (ERA), 389–390
- Estrogens  
  endometrial biopsy and, 387  
  testosterone and, 862  
  thyroxine-binding globulin and, 897  
  total thyroxine (T<sub>4</sub>) and, 895  
  triglycerides and, 908  
  triiodothyronine (T<sub>3</sub>) radioimmunoassay  
    and, 911
- Estrone (E<sub>1</sub>), 386*t*, 387
- ESV (end-systolic volume), 202–203*t*
- Ethanol (EtOH), 391–392  
  alcohol levels, 391–392  
  thromboelastography and, 874–875
- Ethyl alcohol, 391–392  
  alcohol levels, 391–392
- Event marker, in Holter monitoring,  
510–511
- Event recorder, in Holter monitoring,  
510–511
- Evoked brain potentials, 393–396
- Evoked potential studies (EP), 393–396
- Evoked responses, 393–396
- Examination, of bone marrow, 160–164,  
160*t*, 164*t*
- Excretions  
  of bilirubin, 137, 138*f*  
  urine sodium, 837
- Excretory urography (EUG), 764–768
- Exercise stress testing, 210–213
- Expiratory flow rate, peak, 761
- Expiratory flow<sub>25-75</sub>, forced, 761
- Expiratory flow<sub>200-1200</sub>, forced, 761
- Expiratory reserve volume (ERV), 760*f*, 761
- Expiratory volume in 1 second, forced,  
759–760
- Expired carbon monoxide testing, 647
- External sphincter (skeletal muscle) activity,  
during voiding, 353
- Extractable nuclear antigens, antibodies to,  
75–76, 75–76*t*
- Extratesticular lesions, 799
- Exudates  
  peritoneal, 671  
  in pleural fluid analysis, 864
- F**
- F1 + 2 (prothrombin fragment), 877, 878*t*
- Factor I (fibrinogen), 422–423, 423*t*
- Factor V, 397–398
- Factor V Leiden (FVL), 397–398
- Factor VIII, 248
- Factor XI deficiency, activated partial  
  thromboplastin time and, 682
- Familial adenomatous polyposis (FAP), 447
- Familial hypercholesterolemia type B, apo  
  B and, 102
- Familial medullary thyroid carcinoma  
  (FMTC), 449
- Familial thrombosis, plasminogen activator  
  inhibitor 1 and, 700
- Far-advanced cirrhosis, 748*t*
- Fasting blood sugar (FBS), 462–464, 464*t*
- Fasting hypoglycemia, insulin assay in, 539
- Fat, fecal, 402–403, 403*t*  
  total output of, 402
- Fat absorption, 402–403, 403*t*
- Fat retention coefficient, 402–403
- Fatty casts, 936, 943
- Febrile agglutinins, 399, 399*t*
- Febrile antibodies (agglutinins), 399, 399*t*
- Fecal calprotectin, 400–401, 401*t*
- Fecal elastase, 660–661, 661*t*
- Fecal fat, 402–403, 403*t*  
  total output of, 402
- Fecal immunotest (FIT), 844–846
- Fecal incontinence, pelvic floor sphincter  
  electromyography for, 353
- Fecal lactoferrin, 555
- Fecal leukocyte stain, 849
- Fecal leukocytosis, 849
- Fecal occult blood test (FOBT), 844–846
- Female genital tract, pelvic ultrasonography  
  in, 685
- Ferning, of cervical mucus, 822
- Ferritin, 404–405, 405*t*, 546
- Fertility  
  lutinizing hormone in, 593  
  testing, endometrial biopsy for, 387
- Fertilization, sperm antibodies in, 93
- Fetal acid-base status, fetal scalp blood pH  
  and, 417
- Fetal activity determination, 415–416
- Fetal biophysical activity, 406
- Fetal biophysical profile, 406–408
- Fetal blood sampling, fetoscopy for, 420*f*
- Fetal body movements, 407
- Fetal body wall defects, alpha-fetoprotein  
  in, 41
- Fetal breathing movements, 406–407
- Fetal contraction stress test, 409–411
- Fetal distress, 51  
  fetal scalp blood pH and, 417
- Fetal fibronectin (fFN), 412
- Fetal heart rate (FHR), 406  
  reactivity, 406
- Fetal hemoglobin, 413–414
- Fetal hypoxia, 417

## 1036 index

- Fetal infections, in parvovirus B19 antibody, 684
- Fetal lung maturity (FLM), 49  
test in, in L/S ratio, 50
- Fetal maturity, in amniocentesis, 49–50
- Fetal nonstress test, 415–416
- Fetal nuchal translucency (FNT), 612–613, 685
- Fetal oxygen saturation monitoring (FSpO<sub>2</sub>), 417
- Fetal red blood cells, leakage of, 413
- Fetal response, reactive or nonreactive, 415
- Fetal scalp blood pH, 417–418
- Fetal skin biopsies, for primary skin disorders, 419
- Fetal tone, 407
- Fetal-maternal hemorrhage (FMH), 413  
massive, 413  
risk factors for, 413
- a<sub>1</sub>-Fetoprotein, 41–42, 42*t*
- Fetoscopy, 419–421, 420*f*
- Fetus  
direct visualization of, 419, 420*f*  
ultrasound of, 685
- Fiberoptic bronchoscope, 185–186  
flexible, 183
- Fiberoptic duodenoscope, 365, 366*f*
- Fiberoptic endoscope, ERCP and, 364
- Fibrates, triglycerides and, 908
- Fibrin, 195–196*t*
- Fibrin clot formation, 875*f*
- Fibrin degradation products (FDPs), 319–320, 330, 877–878, 878*t*
- Fibrin monomers, 877–878, 878*t*
- Fibrin split products (FSPs), 319–320, 877–878, 878*t*
- Fibrin stabilizing factor, 252*t*
- Fibrinogen, 195–196*t*, 248, 249–250*t*, 252*t*, 422–423
- Fibrinolysin, 699, 699*t*
- Fibrinolysis  
hemostasis and, 251*f*  
plasminogen activator inhibitor 1 and, 700
- Fibrinolytic therapy, 699
- Fibrinopeptide A (FPA), 877, 878*t*
- Fibronectin, fetal, 412
- Fibrosis, cystic, 857
- Fifth disease, 684
- Fine needle aspiration biopsy (FNAB), thyroid, 881–882, 881*t*
- Fine needle aspiration (FNA), 881–882, 881*t*
- First trimester, screening during, for genetic defects, 613
- Fistulas, colon and, barium enema in, 128
- Flexible fiberoptic bronchoscope, 183
- Flow cytometric immunophenotyping, peripheral blood of, 697
- Flow cytometric method, for fetal hemoglobin determination, 413
- Flow cytometry cell surface immunophenotyping, 219–221, 219*t*, 221*t*
- Fluid analysis, 1005–1006  
pancreatic enzymes and, 662  
paracentesis and, 671  
pericardiocentesis and, 689  
semen analysis and, 801  
sexual assault testing of, 809  
Sims-Huhner test and, 822  
sweat electrolytes test and, 857
- Fluorescein, 424
- Fluorescein angiography (FA), 424–425, 425*t*
- Fluorescence in situ hybridization (FISH), 665  
in bladder cancer, 145  
for neoplastic cervical disease, 264
- Fluorescent treponemal antibody test (FTA), 859–860
- Foam stability index test, 50
- Folate, 426–427, 427*t*  
deficiency of, 426  
dietary deficiency of, in homocysteine, 513
- Folic acid, 426–427, 427*t*
- Follicle-stimulating hormone (FSH) assay, 593–594, 593–594*t*
- Forced expiratory flow<sub>25-75</sub> (FEF<sub>25-75</sub>), 761
- Forced expiratory flow<sub>200-1200</sub> (FEF<sub>200-1200</sub>), 761
- Forced expiratory volume in 1 second (FEV<sub>1</sub>), 759–760
- Forced midexpiratory flow, 760
- Forced vital capacity (FVC), 759
- Forensic genetic testing, 445–452, 446*t*
- Four-dimensional computed tomography (4DCT), 281–282
- Fractional excretion of sodium (FE<sub>Na</sub>), 837
- Fractionated metanephrines, 621–622, 622*t*
- Fractionation, of total bilirubin, 138–139
- Fracture assessment, vertebral, 158
- Fragility, erythrocyte, 371–372, 372*t*
- Fragment D-dimer, 319–320
- Francisella tularensis*, in tularemia, 142*t*, 143
- Free cortisol, 290–292, 292*t*
- Free erythrocyte protoporphyrin (FEP), 717
- Free PSA, 743
- Free thyroxine (FT<sub>4</sub>), 894–896. *See also* Thyroxine (T<sub>4</sub>)
- Friedewald formula, 566



Frontal plane view, in 12-lead ECG, 342, 343f  
 Front-line testing measures, 967  
 Fructosamine, 471  
 FTA absorption test (FTA-ABS), 859  
 Fully integrated screen test, 612  
 Functional antithrombin III assay, 97–98, 98t  
 Functional residual capacity (FRC), 760f, 761  
 Functioning nodule, in thyroid scanning, 883  
 Fungal antigen assay, 428–429, 429t  
 Fungal culture, 119, 428–429, 429t  
 Fungal PCR testing, 428–429, 429t  
 Fungal testing, 428–429, 429t  
 Fungi, peritoneal fluid and, 671–672  
 Fungitell, 428–429, 429t  
 Furosemide, free thyroxine ( $T_4$ ) and, 895  
 Furosemide (Lasix) renal scan, 781  
 Fusion CT/positron emission tomography (PET), 273  
 FVL-associated thrombophilia, factor V Leiden for, 397

## G

Gadolinium, 604  
 Gadolinium-containing contrast media, aluminum and, 43  
 Gadoxetate, 601  
 Gail breast cancer risk models, 179  
 GAL-3 (galectin-3), 430, 430t  
 Galactosemia, newborn metabolic screening and, 643  
 Gallbladder  
   nuclear scanning of, 431–432  
   ultrasound of, 1–3  
 Gallium scan, 433–434  
 Gamma activity, decreased, 779  
 Gamma detector/camera, in bone densitometry, 159  
 Gamma-glutamyl transferase (GGT), 435–436, 436t  
 Gamma-glutamyl transpeptidase (GGTP), 435–436, 436t  
 Gangliosidosis, 447  
*Gardnerella vaginalis*, 813t  
 Gas dilution tests, 761–762  
 Gastric atrophy, 688  
 Gastric emptying scan, 437–438, 437t  
 Gastric mucosal biopsy, in *Helicobacter pylori* testing, 482  
 Gastrin, 439–440, 440t  
 Gastrin stimulation tests, 439  
 Gastrin-producing pancreatic tumor, 439  
 Gastroesophageal (GE) reflux  
   esophageal function studies for, 377  
   scan, 441–442, 442t  
 Gastrografin, 922  
   in barium enema, 128  
   for small bowel follow-through, 830  
 Gastrointestinal anthrax, 141  
 Gastrointestinal (GI) bleeding scan, 443–444, 444t  
 Gastrointestinal (GI) scintigraphy, 443–444, 444t  
 Gastrointestinal system, tests for, 991–992  
 Gastroscope, in *Helicobacter pylori* testing, 482  
 Gastroscopy, 382–385  
 Gated blood pool scan, 208  
 Gaucher disease, angiotensin-converting enzyme in, 61  
 G-cell hyperplasia, 439  
 General health, tests for, 1012  
 Genetic aberrations, 51  
 Genetic testing, 445–452, 446t  
   for DNA mismatch repair, 256–257  
   forensic, 445–452, 446t  
 Genital herpes, in sexually transmitted disease testing, 813–814  
 Genital HPV infection, HPV test for, 520  
 Genital HPV strains, 520  
 Genomic prostate assay, 744  
 Genomic testing, 453–457  
   second, 454  
 Genomics, 453–457  
 Genotypic testing, of HCV, 496–497  
 Genotyping, 453–457  
 Gentamicin, 333t  
 German measles test, 793–794, 793t  
 Gestational age, fetal biophysical profile and, 406  
 Gestational diabetes mellitus (GDM), 465  
 GI scintigram, 443  
 GI series, lower, 128–186  
*Giardia* spp. (protozoans), in stool, 847  
 Gigantism, 474  
 Gliadin antibodies, 458–459, 458t  
 Globulin, 746–750, 750t  
   in cerebrospinal fluid, 578  
   ratio, 747  
 Glomerular basement antibody, 79  
 Glomerular disease, 628  
 Glomerular filtration rate (GFR), 301  
   estimated, 303–305, 304–305t  
 Glucagon, 460–461, 461t, 692  
   ERCP and, 366  
 Glucagon provocative test, 693  
 Glucagon stimulation test, 692–693  
 Glucagonoma, 460  
 (1,3)- $\beta$ -D-glucan, 428–429, 429t  
 Glucose  
   anaerobic metabolism of, 551  
   blood, 462–464, 464t

## 1038 index

- Glucose (*Continued*)  
in cerebrospinal fluid, 578  
mean plasma, 472, 472*t*  
peritoneal fluid and, 671–672  
in pleural fluid analysis, 866  
postprandial, 465–466, 466*t*  
urinary, 934, 941
- Glucose screen, 1-hour, 465–466, 466*t*
- Glucose tolerance test (GTT), 467–470, 468*f*  
insulin assay in, 539
- Glucosuria, 934
- Glutamic acid decarboxylase antibody (GAD Ab), 327, 327*t*
- Glutamine, in cerebrospinal fluid, 579
- Glycated albumin, 471
- Glycated proteins, 471
- Glycohemoglobin, 471–473, 472–473*t*
- Glycolated hemoglobin, 471–473, 472–473*t*
- Glycolysis, anaerobic, fetal hypoxia causing, 417
- Glycosylated hemoglobin (GHb), 462–463, 471–473, 472–473*t*
- Glycosylation, 471
- Gonadal tumors, adrenal steroid precursors in, 19
- Gonorrhea, 813*t*, 814  
in sexually transmitted disease testing, 813–814
- Goodpasture antibody, 79
- Goodpasture syndrome, 79
- Gout, uric acid and, 927
- Graft-*versus*-host disease (GVHD), 154
- Gram negative bacteria, 981
- Gram positive bacteria, 981
- Gram stain  
of cerebrospinal fluid, 578  
in pleural fluid analysis, 866  
sputum and, 841–842
- Granular casts, 936, 943
- Granulocyte antibodies, 639
- Granulocytes, 149, 975
- Granulomatous diseases, angiotensin-converting enzyme in, 61
- Granulomatous vascular diseases, antineutrophil cytoplasmic antibody in, 84
- Graves disease, 889
- Group AB blood, 153
- Group O blood, 153
- Growth factors, 314
- Growth hormone (GH), 474–476, 476*t*, 541
- Growth hormone (GH) provocation test, 477–478
- Growth hormone (GH) stimulation test, 474–478, 476*t*
- Growth hormone (GH) suppression testing, 475
- Guaiaac, 844
- Guanethidine  
catecholamines and, 956  
vanillylmandelic acid and, 955
- ## H
- Hageman factor, 249–250*t*, 252*t*
- HAI rubella titer, 794
- Hair, for drug testing, 852
- Haptoglobin, 479–480
- Hashimoto thyroiditis, antithyroid peroxidase antibody in, 100
- hCG beta subunit, 515–517
- HE4, as cancer tumor markers, 195–196*t*
- Health, general, tests for, 1012
- Heart  
computed tomography of, 284–286  
electrocardiography of, 342–346, 344*f*  
electrophysiologic study of, 359–361  
MRI of, 600  
phase-contrast magnetic resonance imaging of, 600  
sonogram of, 339–341
- Heart conduction system, defects in, 359
- Heart failure, 430
- Heart rate  
fetal, 406  
target, 210–211
- Heart scan, 207–209
- Heart valves, artificial, 770
- Heavy dysfunctional uterine bleeding, hysteroscopy in, 534
- Heavy metal (lead) intoxication, 717
- Heinz body preparation, 481
- Helical/spiral CT scan  
of abdomen and pelvis, 272–277, 277*t*  
of brain, 278–280, 280*t*  
of chest, 281–283
- Helicobacter pylori*, 925
- Helicobacter pylori* breath test, 925–926
- Helicobacter pylori* stool antigen, 482–484
- Helicobacter pylori* testing, 482–484
- Hematocrite, 799
- Hematocrit (Hct), 485–487, 486*f*  
hemoglobin and, 488
- Hematologic cancer genomics, 456
- Hematologic diseases, bone marrow examination of, 160
- Hematologic system, tests for, 992–993
- Hematoma, haptoglobin in, 479
- Heme, 717  
bilirubin and, 137
- Hemocult slide test, of stool, 846
- Hemochromatosis, 445–452, 446*t*, 547  
genetic testing of, 449

- Hemochromatosis-associated HFE genes, 449
- Hemoglobin A<sub>1c</sub> (HbA<sub>1c</sub>), 471, 472*t*
- Hemoglobin C (Hgb C), 490, 491*t*
- Hemoglobin E (Hgb E), 490, 491*t*
- Hemoglobin F (Hgb F), 490, 491*t*
- Hemoglobin H (Hgb H), 491*t*
- Hemoglobin (Hb, Hgb), 488–489, 772
- Apt test in, 105
- bilirubin and, 137
- content of, in hemoglobinopathies, 491*t*
- electrophoresis, 490–492
- fetal, 413–414
- glycolated, 471–473, 472–473*t*
- Hemoglobin (Hgb) S, 490, 491*t*, 820
- test, 820–821
- Hemoglobin M disease, 623
- Hemoglobinopathies, hemoglobin in, 488
- contents of, 491*t*
- electrophoresis of, 490
- Hemolysis, 287
- in ammonia, 47
- haptoglobin in, 479
- Heinz bodies in, 481
- osmotic fragility in, 371
- profile, 1012
- Hemolytic anemias
- haptoglobin in, 479
- Heinz bodies in, 481
- Hemolytic disease, Rh typing in, 153
- Hemorrhage, fetal-maternal, 413
- Hemorrhagic fever, bioterrorism infectious agents testing for, 141, 142*t*
- Hemosiderin, 546
- Hemosiderosis, 547
- Hemostasis, 874
- coagulation factors concentration and, 248
- fibrinolysis and, 251*f*
- Heparin
- in activated clotting time, 10
- in antithrombin activity, 97
- cofactor, 97–98, 98*t*
- free thyroxine (T<sub>4</sub>) and, 895
- platelet antibody detection and, 704
- therapeutic actions, 681
- thromboclastography and, 874–875
- thyroid-stimulating hormone and, 887
- WBC levels and, 976
- Heparin-induced thrombocytopenia antibodies (HITAs), 704
- Heparin-induced thrombocytopenia (HIT), 704
- Hepatic cellular injury, aldolase in, 23
- Hepatic encephalopathy, ammonia in, 47
- Hepatic function panel, 1012
- Hepatic ultrasound, 1–2
- Hepatitis, 493, 813*t*
- acute, tests for, 1012–1013
- testing of, 495*t*
- viral, ALT/AST ratio in, 21
- Hepatitis A virus (HAV), 493
- RNA, 495*t*
- Hepatitis B core antibody (HBcAb), 494, 495*t*
- Hepatitis B core antigen (HBcAg), 494, 495*t*
- Hepatitis B e-antibody (HBcAb), 494–496
- Hepatitis B e-antigen (HBcAg), 494–496, 495*t*
- Hepatitis B surface antibody (HBsAb), 494
- Hepatitis B surface antigen (HBsAg), 494–496, 495*t*
- Hepatitis B virus (HBV), 493–494
- DNA, 495*t*, 496
- viral load, 496
- Hepatitis C, anti-LKM-1 antibodies in, 81
- Hepatitis C antibody (HCV-Ab), 495*t*
- Hepatitis C virus (HCV), 496
- genotypic testing of, 496–497
- RNA, 495*t*
- testing, 496
- Hepatitis D antibody (HDV-Ab), 495*t*
- Hepatitis D antigen (HDV-Ag), 495*t*, 497
- Hepatitis D virus (HDV), 497
- Hepatitis E virus (HEV), 497
- Hepatitis virus studies, 493–497
- Hepatobiliary disease, prealbumin and, 728–729
- Hepatobiliary iminodiacetic acid analogue (HIDA) scanning, 431–432
- Hepatobiliary scintigraphy, 431–432
- Hepatobiliary system, tests for, 993
- Hepatocellular disease
- ALT/AST ratio in, 21
- ammonia in, 47
- aspartate aminotransferase in, 125
- Hepatocellular dysfunction, in jaundice, 137–138
- Hepatocellular liver disease, 753
- HER2, 588
- testing for, 588–589
- Hereditary AT-III deficiency, 97
- Hereditary hemochromatosis (HH), 449
- Hereditary nonpolyposis colorectal cancer (HNPCC), 256
- Hereditary nonpolyposis colorectal cancer (HNPCC) syndrome, 447
- Hereditary spherocytosis, erythrocyte fragility for, 371
- Heroin, total thyroxine (T<sub>4</sub>) and, 895
- Herpes, TORCH test for, 899
- Herpes genitalis, 498–499, 813*t*
- Herpes simplex, 498–499
- Herpes simplex virus, in sexually transmitted disease testing, 813–814

## 1040 index

- Herpes simplex virus type 1 (HSV 1), 498–499
- Herpes simplex virus type 2 (HSV 2), 498–499
- Herpes virus-3, 958–959
- Herpesvirus type 1, 498–499
- Herpesvirus type 2, 498–499
- Heterophil antibody test, 630
- Heterophile antibodies  
for EBV infections, 368, 370  
tuberculosis testing and, 920
- Heterozygous individuals, for FVL, 397
- Hexosaminidase A (Hex A), 500–501
- Hexosaminidase B (Hex B), 500–501
- Hexosaminidase (HEX), 447, 500–501
- HFE genes, hemochromatosis-associated, 449
- Hgb A<sub>1</sub>, 490, 491*t*
- Hgb A<sub>2</sub>, 490, 491*t*
- High-density lipoproteins (HDLs), 235, 565–568, 567–568*t*  
apolipoproteins in, 101–102  
subclasses of, 565–566
- High-grade cervical intraepithelial lesions, in HPV, 520
- High-sensitivity C-reactive protein (hs-CRP), 295–296, 296*t*
- High-sensitivity troponin (hsTnT, hs-cTn), 912–914. *See also* Troponins
- Histocompatibility leukocyte A antigen, 518–519
- Histone antibodies, 68
- HIV antigen/antibody (Ag/Ab)  
combination assays, 508
- HIV-RNA viral test, 506–509
- Holter monitoring, 510–512, 510*f*
- Home sleep studies, 827–828
- Homocysteine (Hcy), 513–514
- Homocystinuria, newborn metabolic screening and, 644
- Homologous blood, for donation, 154
- Homovanillic acid (HVA), 954–957
- Homozygous individuals, for FVL, 397
- Horizontal plane view, in 12-lead ECG, 342, 343*f*
- Hormone assays, for urinary pregnanediol, 731
- Hormone-producing tumors, 652
- Hot spots, 165
- Hot water, caloric study and, 192
- 1-Hour glucose screen, 465–466, 466*t*
- 2-Hour postprandial glucose (2-Hour PPG), 465–466, 466*t*
- 24-hour test, in adrenocorticotropic hormone stimulation test with cosyntropin, 14
- Howell-Jolly bodies, in red blood cells, microscopic examination of, 150*b*
- Hughes syndrome, anticardiolipin antibodies in, 65
- Human blood, group of, 152
- Human calcitonin (HCT), 187–188, 188*t*
- Human chorionic gonadotropin (hCG, HCG), 195–196*t*, 515–517
- Human chorionic somatomammotropin (HCS), 523–524
- Human diploid cell rabies vaccine (HDCV), 769
- Human epidermal growth factor receptor 2 (HER-2), protein, 176–178
- Human growth hormone (HGH), 474–476, 476*t*
- Human immunodeficiency virus (HIV), 813*t*  
antibody test, 506–509  
Centers for Disease Control and Prevention, 507*t*  
drug resistance testing, 502  
drug sensitivity testing, 502  
genotype, 502  
HTLV and, 525  
risk factors for, 507*b*  
RNA quantification, 503–505  
serologic and virologic testing, 506–509, 507*t*  
sexually transmitted disease testing, 813–814  
tests for, 1013  
viral load, 503–505
- Human immunodeficiency virus type 1 (HIV 1), 506
- Human immunodeficiency virus type 2 (HIV 2), 506
- Human interferon inducible protein 10, 314
- Human leukocyte A antigen, 518–519
- Human lymphocyte antigen B27 (HLA-B27 antigen), 518–519
- Human lymphocyte antigen (HLA), 704
- Human papillomavirus (HPV), 520–522, 521*t*  
in sexually transmitted disease testing, 813–814
- Human placental lactogen (hPL), 523–524
- Human T-cell lymphotropic virus (HTLV) I/II antibody, 525
- Huntington disease, PET scan and, 719–720
- Hyaline casts, 936, 943
- Hyaluronan, 803
- Hyaluronan binding assay (HBA), 803
- Hydrocele, 799
- Hydrocortisone, 290–292, 292*t*
- Hydrogen breath test, 557, 559*t*
- Hydrogen (H<sub>2</sub>), in SIBO test, 833
- Hydronephrosis, 776
- 3-Hydroxy cotinine, 646–648

- 25-Hydroxy vitamin D<sub>2</sub>, 971–973  
 25-Hydroxy vitamin D<sub>3</sub>, 971–973  
 17-Hydroxycorticosteroids (17-OCHS),  
 16, 526–527  
 5-Hydroxyindoleacetic acid (5-HIAA),  
 528–529  
 21-Hydroxylase antibodies, 530  
 17-Hydroxypregnenolone, 18–19, 18–19<sup>t</sup>  
 17-Hydroxyprogesterone, 18–19, 18–19<sup>t</sup>  
 5-Hydroxytryptamine (5-HT), 807–808,  
 808<sup>t</sup>  
 Hyperaldosteronism, 783–784  
 aldosterone in, 25  
 Hyperamylasemia, amylase in, 56  
 Hyperbilirubinemia, 935  
 unconjugated, 138–139  
 Hypercalcemia, 189–190, 675, 677  
 Hyperchloremia, 233  
 Hyperchromasia, in red blood cells,  
 microscopic examination of, 150<sup>b</sup>  
 Hyperchromic anemia, 774  
 Hypercoagulable disorders, antithrombin  
 activity in, 97  
 Hypercoagulable state, 874, 876  
 Hypercortisol, 324  
 Hyperfunctioning adrenal tumors, Cushing  
 syndrome and, 13  
 Hyperglycemia, insulin assay in, 539  
 Hyperkalemia, 724  
 Hypermagnesemia, 597  
 Hypermnatremia, 836<sup>t</sup>  
 Hyperparathyroidism  
 parathyroid hormone, 675  
 parathyroid scan, 677  
 Hyperphosphatemia, 694  
 Hyperplasia, G-cell, 439  
 Hypersensitivity reaction, in blood  
 transfusion, 153  
 Hyperstimulation, uterine, 409  
 Hypertension  
 portal, 574, 770  
 renovascular, 783, 784<sup>f</sup>  
 renin assay, 783–784  
 tests for, 1013  
 total blood volume for, 900  
 Hyperthyroidism, 892, 910  
 Hypertonia, human placental lactogen  
 in, 523  
 Hypertrophic cardiomyopathy  
 (HCM), 450  
 Hyperuricemia, 927  
 Hyperventilation, 349  
 Hypocalcemia, 190  
 Hypochloremia, 233  
 Hypochromic anemia, 774  
 Hypocoagulable state, 874, 876  
 Hypogammaglobulinemia, 748<sup>t</sup>  
 Hypokalemia, 724  
 Hyponatremia, 836<sup>t</sup>  
 aldosterone in, 25  
 water load test in, 72  
 Hypoparathyroidism, 675  
 Hypophosphatemia, 694  
 Hypopituitarism, adrenocorticotrophic  
 hormone stimulation test with  
 cosyntropin in, 13  
 Hypoproteinemia  
 anion gap in, 63  
 transferrin in, 547  
 Hypothyroidism, congenital, newborn  
 metabolic screening and, 643  
 Hypovolemia, 672  
 Hypoxia, 551  
 fetal, 417  
 Hysteroqram, 531–532  
 Hysterosalpingography, 531–532  
 Hysteroscopy, 533–535, 533<sup>f</sup>
- I**  
 Idiopathic thrombocytopenia purpura  
 (ITP), 704  
 IgA anti-*H. pylori* antibody, 483  
 IgA antisperm antibodies, 93  
 IgE antibody test, 31–32, 32<sup>t</sup>  
 IgG antibody  
 anti-*H. pylori* antibody, 483  
 for HAV, 493  
 in *Helicobacter pylori* testing, 482<sup>t</sup>  
 IgG antibody testing, 31–32  
 IgG antisperm antibodies, 93  
 IgG (immunoglobulin G)-specific  
 antibodies, 684  
 IgM (immunoglobulin M) antirubella  
 antibodies, 794  
 IgM (immunoglobulin M) titer, 794  
 Ileocolic intussusception, nonstrangulated,  
 barium enema in, 128  
 Iminodiacetic acid analogs (IDAs), 431  
 Imipramine, vanillylmandelic acid and,  
 955  
 Immunochemical fecal occult blood test  
 (iFOBT), 844  
 Immunofixation, 134  
 in immunoglobulin quantification, 536  
 Immunofixation electrophoresis (IFE),  
 746–750, 750<sup>t</sup>  
 Immunofluorescence antibody (IFA), 824  
 Immunoglobulin A (IgA), 458, 536  
 Immunoglobulin D (IgD), 536  
 Immunoglobulin E (IgE), 536  
 Immunoglobulin G (IgG), 536  
 Immunoglobulin light chains, 134  
 Immunoglobulin M (IgM), 536  
 negative for, 684  
 reactive, 790

## 1042 index

- Immunoglobulin M (IgM) antibodies  
anti-*H. pylori* antibody, 483  
of HAV, 493  
in *Helicobacter pylori* testing, 482*t*  
Zika virus and, 985
- Immunoglobulin quantification, 536–538, 537–538*t*
- Immunoglobulins, 195–196*t*
- Immunologic antithrombin III, 97–98, 98*t*
- Immunologic system, tests for, 994–995
- Immunologic tests, 368
- Immunophenotyping, cell surface, 219–221, 219*t*, 221*t*
- Immunosuppressed patients, EBV infection and, 370
- Impedance, chest, 826
- Implantable cardiac defibrillators (ICDs), MTWA and, 344–345
- Implantable loop recorders (ILRs), 511
- In vitro test, for specific IgE, 31–32
- Inactivate prothrombin, 681
- Inactive (dormant) sarcoidosis, angiotensin-converting enzyme in, 61
- Inborn errors of metabolism, 642
- Indirect antiglobulin test (IAT), 289
- Indirect bilirubin, 137
- Indirect Coombs test, 289
- Industrial testing, 852
- Indwelling urinary catheter, 948
- Infections  
amniotic fluid for, assessment of, 51  
fungal, 428  
parvovirus B19 antibody  
acute, 684  
fetal, 684  
in prosthetic joints, 39  
timing of  
EBV antibodies and, 369*t*  
serologic studies and, 369*t*
- Infectious agents testing, bioterrorism, 141–144, 142*t*
- Infectious diseases, febrile antibodies for, 399
- Infectious hepatitis, 493
- Infectious mononucleosis, 630
- Infertility  
male, semen analysis and, 802  
screen, 93–94  
Sims-Huhner test and, 822  
in women, 593
- Inflammatory arthritis, viscosity in, 118
- Inflammatory autoimmune diseases, erythrocyte sedimentation rate for, 373
- Inflammatory bowel disease, barium enema in, 128
- Inflammatory diseases, haptoglobin in, 479
- Inflammatory intestinal disorders, lactoferrin and, 555
- Inflammatory scan, 979–980
- Infusion test, calcium, 187
- Ingested fat, percentage of, 402–403
- Inhalant allergens, allergy skin testing in, 34
- Inhalation tests, 763
- Inherited AAT deficiency, 37
- Inhibin A, 195–196*t*, 613
- Inorganic phosphate, 694
- Inspiratory capacity (IC), 760*f*, 761
- Inspiratory flow rate, peak, 761
- Inspiratory reserve volume (IRV), 760*f*, 761
- Insulin, catecholamines and, 956
- Insulin assay, 539–540, 540*t*
- Insulin autoantibody (IAA), 327, 327*t*
- Insulin C-peptide, 293–294, 294*t*
- Insulin tolerance test (ITT), 477–478
- Insulin-dependent diabetes mellitus (IDDM), 327
- Insulin-like growth factor (IGF), 541–542, 541–542*t*
- Insulin-like growth factor binding proteins (IGF BP), 541–542, 541–542*t*
- Insulin-like growth factor (IGF-1), 474, 541–542
- Insulinoma, insulin assay in, 539
- Interferon, 314
- Interferon gamma release assay (IGRA), 919–921. *See also* Tuberculosis (TB), testing
- Interferon inducible protein 10, human, 314
- Interleukins, 314–315
- International normalized ratio (INR), 753–756, 754*t*, 756*t*
- Intestinal infarction, elevated lipase level and, 562
- Intracoronary stents, 203
- Intradermal test, in allergy skin testing, 33, 35
- Intrauterine device (IUD), 686
- Intravascular catheter, percutaneous IVUS imaging in, 543
- Intravascular hemolysis, osmotic fragility and, 371
- Intravascular thrombosis, protein C, protein S, 751
- Intravascular ultrasound (IVUS), 543–544
- Intravenous pyelography (IVP), 764–768
- Intravenous urography (IUG), 764–768
- Intrinsic clotting pathway, 681
- Intrinsic factor antibody (IF ab), 545, 545*t*
- Intrinsic factor (IF), 969
- Intussusception, nonstrangulated ileocolic, barium enema in, 128
- Iodinated x-ray contrast dye, 364
- Iodine-containing contrast media, aluminum and, 43

Iodine-123 capsule, in thyroid scanning, 883

Iontophoretic sweat test, 857–858

Iron, ferritin storage of, 404

Iron deficiency anemia, 404, 547  
zinc protoporphyrin for, 987

Iron level, 546–549, 549*t*  
abnormal findings in, 549  
contraindications for, 548  
interfering factors in, 548  
normal findings in, 546  
procedure and patient care for, 548  
test explanation and related physiology of, 546–548

Iron poisoning, acute, 547

Iron preparations (Ferotran, Imferon), urine color and, 932*t*

Irritable bowel syndrome (IBS), 316  
lactoferrin and, 555

Ischemia, ischemia-modified albumin in, 550

Ischemia-modified albumin (IMA), 550, 550*t*

Ischemic stroke, PLAC test for, 564

Islet cell antibody (ICA), 327, 327*t*

Isoenzymes, 297, 553  
of alkaline phosphatase, 29

Isonitrile scan, 207–209

*Ixodes pacificus*, 595

*Ixodes scapularis*, 595

## J

Jaundice, 137

endoscopic retrograde  
cholangiopancreatography  
for, 364

Joint, tests for, 1011

## K

Kappa light chains, free, 134–135, 135*t*

Karyotype, 244–245, 245*t*

Karyotyping, 241, 244

Ketoacids, 927

Ketoconazole, testosterone and, 862

Ketones, urinary, 935, 938–939, 942

Ketonuria, 935

17-Ketosteroids (17-KS), 862

Ki67 protein, 176–178

Kidney(s)

antidiuretic hormone in, 71  
arteriography, 117  
biopsy of, 776–778, 777*f*  
erythropoietin produced by, 375  
normal, osmolality and, 656  
nuclear imaging of, 779–782  
scan, 779–782  
ultrasound of, 1

Kidney disease

and B<sub>2</sub>m, 628

microalbumin and, 626

Kidney stones, 928

Kleihauer-Betke test, 413–414

KRAS, 588

as cancer tumor markers, 195–196*t*

testing for, 588–589

KUB (kidney, ureter, bladder) film, x-ray of, 650

## L

Lactase deficiency, 557

Lactate, 551–552

Lactate dehydrogenase (LDH), 195–196*t*, 350

Lactic acid, 551–596, 552*t*

in cerebrospinal fluid, 579

Lactic acidosis (LA), 551, 552*t*

brain CT scan in, 279

chest CT scan in, 282

metformin and, 115

type I, 551

type II, 551

type III, 551

Lactic dehydrogenase (LDH), 553–554, 554*t*

in cerebrospinal fluid, 579

isoenzyme, 553

peritoneal fluid and, 671–672

in pleural fluid analysis, 866

Lactoferrin, 555–556

fecal, 555

white blood cells, 555

Lactose breath test, 557

Lactose tolerance test, 557–559, 559*t*

Lactulose breath test, 833–834

Lambda light chains, 134–135, 135*t*

Lamellar body count, 50

Large loop excision of transformation zone (LLETZ), 224

Large-volume blood transfusions, fibrinogen and, 422

Laron dwarfism, IGF BPs in, 542

Laryngoscopy, 184

Late-onset Tay-Sachs disease, 447

Late-phase reaction, in allergy skin testing, 33

Lateral view, in chest x-ray, 227

Lavage, breast ductal, 179–180

Lead, 620, 620*t*

12-Lead ECG, 342, 343*f*

Lead intoxication, porphyrins and porphobilinogens, 717

Lead poisoning, zinc protoporphyrin and, 987

Lecithin, 49–50

## 1044 index

- Lecithin/sphingomyelin (L/S ratio), in amniocentesis, 49–50
- Left ventricular pressure, end-diastolic, 202–203*t*
- Legionella pneumophila*, 560
- Legionnaires disease, 560
- antibody test, 560, 560*t*
  - diagnosis of, 560
- Leiden factor V, 397–398
- Lesions, extratesticular, 799
- Leucine aminopeptidase (LAP), 561, 561*t*
- Leukemia, 770
- bone marrow biopsy and, 161
- Leukoagglutinin, 639
- Leukocyte count, 974–978, 974*t*, 977–978*t*
- abnormal findings in, 977–978, 977*t*
- Leukocyte esterase, 935, 938, 942
- Leukocytes, 849
- peripheral blood smear in, 149
  - stool for, 849
  - in synovial fluid, 119
- Leukocytosis, 974
- fecal, 849
- Leukopenia, 974
- Levodopa
- urine color and, 932*t*
  - vanillylmandelic acid and, 955
- Lidocaine, 333*t*
- Light chains, 134
- Lipase, 562–563, 662–664, 664*t*
- Lipid fractionation, 565–568
- Lipid panel, 1013
- Lipid profile, 565–568
- testing, 235
- Lipocalin-2, 640–641, 641*t*
- Lipoprotein (a), 101–103
- Lipoprotein electrophoresis, 565–568
- Lipoprotein phenotyping, 565–568
- Lipoprotein-associated phospholipase A<sub>2</sub>, 564, 564*t*
- Lipoproteins, 565–568
- high-density, 235, 565–568, 567–568*t*
    - apolipoproteins in, 101–102
    - subclasses of, 565–566  - low-density, 235, 565–568, 567–568*t*
    - apolipoproteins in, 101–102
    - subclasses of, 566  - very-low-density, 101–102, 565–568
- Liquid biopsy, 569–570
- Liquid-based cervical cytology (LBCC), 667–670, 670*t*
- Lithium, 333*t*
- thyroid-stimulating hormone and, 887
  - total thyroxine (T<sub>4</sub>) and, 895
  - vanillylmandelic acid and, 955
- Liver
- 5' nucleotidase to, 649
  - bile in, 137
- Liver (*Continued*)
- biopsy of, 571–573, 571*f*
    - contraindications to, 572
    - indications for, 571
    - potential complications of, 572
    - in primary biliary cirrhosis, 82  - disorders of, LAP in, 561
  - fibrinogen produced by, 422
  - gamma-glutamyl transpeptidase in, 435
  - MRI of, 601
  - nuclear scan of, 574
  - scanning of, 574–575
    - contraindications to, 574  - ultrasound of, 1–2, 4
  - uroporphyrinogen-1-synthase in, 953
- Liver disease, 709
- alkaline phosphatase in, 29
  - primary, haptoglobin in, 479
- Liver parenchyma, in alanine aminotransferase, 21
- LOC387715, in age-related macular degeneration, 20
- Long bone, x-ray in, 172
- Long QT syndrome (LQTS), 450
- Long-acting thyroid stimulator (LATS), 889–890
- Loop electrosurgical excision procedure (LEEP), 224–226
- Lordotic view, in chest x-ray, 227
- Low bone mass, 157
- Low molecular weight heparins (LMWHs), plasma anti-Xa assay in, 77
- Low plasma levels, plasminogen activator inhibitor 1 and, 700
- Low-density lipoproteins (LDLs), 235, 565–568, 567–568*t*
- apolipoproteins in, 101–102
  - subclasses of, 566
- Lower extremities
- arteriography of, 114
  - venography of, 963–964
- Lower GI series, 128–130, 130*t*
- Lp-PLA<sub>2</sub> (lipoprotein-associated phospholipase A<sub>2</sub>), 102, 564, 564*t*
- Lumbar puncture (LP), 576–583, 581*f*
- for cerebrospinal fluid, 576–583
  - myelography and, 632
- Lumbar spine, MRI of, 601
- Lumbar x-ray studies, 839–840
- Lung
- biopsy of, 584–587
    - closed technique of, 584
    - in Goodpasture syndrome, 79
    - open method of, 584, 586  - percutaneous, 843
  - percutaneous needle biopsy as, 586



- Lung (*Continued*)  
 thoracoscopic, 586  
 transbronchial, 585, 585*f*  
 transbronchial brushing, 586  
 transbronchial needle aspiration, 585, 585*f*  
 diffusing capacity of, 763  
 nuclear scan of, 590  
 transfusion-related injury, 639
- Lung cancer  
 genomic testing of, 455  
 molecular testing for, 588–589  
 neuron-specific enolase in, 638
- Lung scan, 590–592
- Lung volumes/capacities, 760*f*
- Lupoid CAH, 92
- Lupus anticoagulant, 65–66, 66*t*
- Lupus erythematosus (LE), ribosome P antibodies for, 792
- Luteinizing hormone, 593–594, 593–594*t*
- Lutropin, 593–594, 593–594*t*
- Lyme disease test, 595–596
- Lymphocyte immunophenotyping, 219–221, 219*t*, 221*t*
- Lymphocyte precursors, bone marrow biopsy in, 162
- Lymphocytes, 975
- Lymphocytopenia, 977–978*t*
- Lymphocytosis, 977–978*t*
- Lymphogranuloma venereum, 813*t*
- Lymphoma cancer genomics, 456
- Lymphoscintigraphy, 805–806
- Lynch syndrome, 447
- Lysosomal storage disease (LSD), 7, 500
- M**
- Macrocytes, in red blood cells, microscopic examination of, 150*b*
- Macrocytic cell, 772
- Magnesium, 597–598, 598*t*  
 deficiency of, 597  
 with potassium, and calcium, 597
- Magnetic resonance angiography (MRA), 600
- Magnetic resonance  
 cholangiopancreatography (MRCP), 364, 601  
 indications for, 601
- Magnetic resonance enterography (MRE), 601
- Magnetic resonance imaging (MRI), 599–605  
 of brain and meninges, 599  
 of breast, 600  
 of cervical and lumbar spine, 601
- Magnetic resonance imaging (MRI) (*Continued*)  
 of heart, 600  
 of liver, 601  
 myelogram, 601  
 upright, 601
- Magnetic resonance spectroscopy (MRS), 599
- Magnetic resonance venography (MRV), 601–602
- Magnetoencephalography (MEG), 347
- Magnified spot views, 610
- Malabsorption, D-xylose absorption test and, 983
- Malignant ventricular arrhythmias, signal-averaged ECG and, 344
- Malnutrition, prealbumin, 729*t*
- Mammary ductoscopy, 336–338, 337*f*, 338*t*
- Mammography, 606–611  
 digital, 606  
 PET, 721  
 screening guidelines, 608–609*t*  
 three-dimensional, 606–607  
 ultrasound, 181–182, 182*t*
- Manometric studies, 364, 377
- Manometric tests, 1006  
 arterial plethysmography and, 713
- Manometry, esophageal, 377–381, 380*f*
- Mantoux test, 915–916
- Maple syrup urine disease (MSUD), newborn metabolic screening and, 643–644
- Marijuana, 853*t*
- Marrow red blood cell (RBC) precursor, bone marrow biopsy in, 161
- Massive fetal-maternal hemorrhage, 413  
 risk factors for, 413
- Mast cell disease, systemic, 11 beta-prostaglandin F(2) alpha in, 136
- Mast cells, 975
- Mastocytosis, systemic, 11 beta-prostaglandin F(2) alpha in, 136
- Masturbation, 804
- Maternal DNA testing, cell-free, 217–218
- Maternal quadruple screen, 612–614
- Maternal screen testing, 612–614
- Maternal serum, PAPP-A in, 730
- Maternal triple screen, 612–614
- Maternal-fetal oxygen transfer, 409
- Maternal-fetal platelet antigen incompatibility, 704
- Maximal breathing capacity, 760
- Maximal midexpiratory flow (MMEF), 760
- Maximal volume ventilation (MVV), 760
- MDMA (ecstasy), 853*t*

## 1046 index

- Mean corpuscular hemoglobin concentration (MCHC), 774
- Mean corpuscular hemoglobin (MHC), 773
- Mean corpuscular volume (MCV), 712, 772
- Mean plasma glucose (MPG) level, 472, 472*t*
- Mean platelet volume (MPV), 712, 712*t*
- Measles rubeola antibody, 615
- Measles virus, 615  
testing for, 615
- Measured osmolality, 654–655
- Measurement, symbols and units of, 1015–1016, 1015–1016*t*
- Meckel diverticulum nuclear scan, 616–617
- Mediastinal lymph nodes, mediastinoscopy for, 618
- Mediastinoscopy, 618–619
- Megakaryocytes, bone marrow biopsy in, 161–162
- Megaloblastic anemias, 772  
haptoglobin in, 479  
homocysteine in, 513
- MEK1, 588  
testing for, 588–589
- Melanoma, 445–452, 446*t*, 805  
genetic testing of, 448–449
- Membrane instability, 371
- Membranous nephropathy (MN), 696
- Meninges, MRI of, 599
- Meningitis, lumbar puncture and cerebrospinal fluid examination, 578, 581*f*
- Mercury, 620, 620*t*
- MET, 588  
testing for, 588–589
- Metabolic acidosis, 108–109*t*
- Metabolic alkalosis, 108–109*t*
- Metabolic profile, basic, 133
- Metabolic screening, newborn, 642–645
- Metabolism, of bilirubin, 137, 138*f*
- Metabolites, nicotine and, 646–648, 646*t*
- Metal testing, 620, 620*t*
- Metanephrine, 954–957  
plasma free, 621–622, 622*t*  
interpretation of, 621  
testing of, 692
- Metastatic neuroendocrine tumors, octreotide scan of, 652
- Metastatic process, in bone scan, 166
- Metformin, 279, 282
- Methodone, 853*t*  
free thyroxine ( $T_4$ ) and, 895  
thyroxine-binding globulin and, 897  
total thyroxine ( $T_4$ ) and, 895  
triiodothyronine ( $T_3$ ) radioimmunoassay and, 911
- Methane ( $CH_4$ ), in SIBO test, 833
- Methemoglobinemia, acquired, 623
- Methenamine, catecholamines and, 956
- Methionine  
loading, 514  
metabolism, inborn errors of, homocysteine in, 513
- Methotrexate, 333*t*
- 3-Methoxy-4-hydroxymandelic acid, 955
- Methylated septin 9 DNA, 625, 625*t*
- Methylation markers, for neoplastic cervical disease, 264
- Methyldopa, catecholamine and, 956
- Methylmalonic acid (MMA), 969–970
- Metronidazole (Flagyl), urine color and, 932*t*
- Metyrapone test, 16–17
- Micral Urine Test Strip, 626
- Microalbumin (MA), 626–627, 627*t*
- Microcytes, in red blood cells, microscopic examination of, 150*b*
- Microcytic cell, 772
- Microglobulin, 628–629, 629*t*
- Microhemagglutination assay (MHA-TP), 859
- Microsatellite instability (MSI) gene, 455
- Microsatellite instability (MSI) testing, 256–257
- Microscopic examination, 1006–1007  
in bone marrow, 161  
pancreatobiliary FISH testing and, 665  
Parkinson disease testing and, 679–680  
of red blood cells, 150*b*  
sexually transmitted disease testing and, 813  
skin biopsy and, 824  
of tissue, pleural biopsy and, 715
- Microsomal antibodies, 100
- Microvolt T-wave alternans (MTWA), 344
- Midstream collection, of urine  
for urinalysis, 939  
for urine culture, 948
- Minimal hemostatic levels, 248
- 40-Minute test period, 415
- Minute volume/ventilation (MV), 761
- Mi-Prostate Score (MiPS), 743–745, 745*t*
- Miscellaneous studies, 995
- Mismatch repair (MMR) gene testing, 455
- Mitotic spindle, centromere in, 67
- Mixed connective tissue disease (MCTD), antiextractable nuclear antigens in, 75
- MM3/MM1 ratio, 299
- M-mode echocardiography, 339
- Modified ACT test, 9–10
- Molecular PCR tissue, for herpes simplex, 499
- Molecular testing  
for lung cancer, 588–589  
Zika virus and, 985

- Monoamine oxidase inhibitors,  
  vanillylmandelic acid, 955
- Monoclonal gammopathy, 748*t*  
  immunoglobulins in, 537
- Monoclonal light chains, 134
- Monocytes, 975
- Monocytopenia, 977–978*t*
- Monocytosis, 977–978*t*
- Mononuclear heterophil test, 630
- Mononucleosis, Epstein-Barr virus testing  
  (EBV) antibody titer for, 368
- Mononucleosis rapid test, 630
- Monospot test, 370, 630
- Mortality, preeclampsia and, 698
- mSEPT9, 625, 625*t*
- M-spike, urine, monitoring, 134
- Mucin clot test, 118–119
- Mucin glycoprotein, 145–146
- Multicenter AIDS Cohort Study, 220
- Multigated acquisition (MUGA) scan,  
  207–209
- Multipanel drug screen, 853*t*
- Multiple electrode catheters, in  
  electrophysiologic study, 359
- Multiple myeloma, Bence Jones protein  
  in, 134
- Multiple sclerosis antibody panel, 80, 80*t*
- Multiple sleep latency test (MSLT), 827,  
  829
- Multiple wake test (MWT), 827, 829
- Multiples of median, 613
- Muscular abnormalities, barium swallow  
  in, 131
- Muscular injury, aldolase in, 23
- Mutation analysis, 397–398
- Mutations, of oncogenes, 588
- Myasthenia gravis (MG), acetylcholine  
  receptor antibody panel in, 5
- Mycobacterium tuberculosis*. *See* Tuberculosis
- Mycoplasma pneumoniae  
  antibodies, IgG and IgM, 631  
  cold agglutinins and, 253
- Myelogram, MRI, 601
- Myelography, 632–634
- Myelomeningocele, alpha-fetoprotein in, 41
- Myeloperoxidase (MPO) autoantigen, 85
- Myeloperoxidase-ANCA (MPO-ANCA), 84
- Myeloproliferative disorders, 702
- Myocardial infarction  
  acute, gamma-glutamyl transpeptidase  
    in, 435  
  creatinine kinase and, 297, 298*f*  
  troponins for, 912–913
- Myocardial injury  
  antimyocardial antibody in, 83  
  lactic dehydrogenase levels in, 553
- Myocardial nuclear stress testing, 207
- Myocardial perfusion imaging, 207–209
- Myocardial perfusion scan, 207–209
- Myocardial scan, 207–209
- Myoglobin, 635, 635*t*
- Myxedema, SIADH and, 72
- ## N
- NAAT (nucleic acid amplification for TB),  
  919–921. *See also* Tuberculosis (TB),  
  testing
- Nasal culture, 873
- Nasopharyngeal carcinoma  
  EA-D and, 368  
  EBV infection and, 368
- Nasopharyngeal culture, 873
- Nasopharyngeal swabs, for SARS viral  
  testing, 797
- Nasopharyngeal wash or aspirates, 797
- Native double-stranded DNA, 74, 74*t*
- Natriuretic hormone, 835
- Natriuretic peptides (NPs), 636–637, 637*t*
- Needle aspiration, transbronchial, 585, 585*f*
- Negative PET, 720
- Negative test, 409
- Neisseria gonorrhoeae*, 814  
  in sexually transmitted disease testing,  
    813–814
- Neonatal morbidity, preeclampsia and, 698
- Neonatal thrombocytopenia, 704
- Neoplastic diseases, associated with febrile  
  agglutinins, 399
- Nephrogenic DI, antidiuretic hormone in,  
  71–72
- Nephrotic syndrome, 748*t*
- Nephrotomography, CT, 274
- Nerve conduction studies (NCS), 355–356
- Nervous system, tests for, 996
- Neural crest tumors, vanillylmandelic acid,  
  955
- Neural tube defect (NTD), 239  
  alpha-fetoprotein in, 41
- Neuroendocrine nuclear scan, 652–653
- Neuroendocrine tumors, 638*t*, 652
- Neurogenic DI, antidiuretic hormone in, 71
- Neurologic diseases, aldolase in, 23
- Neurology, PET scan and, 719–720
- Neuron-specific enolase (NSE), 195–196*t*,  
  638, 638*t*
- Neurosarcoidosis, angiotensin-converting  
  enzyme in, 61
- Neutropenia, 977–978*t*
- Neutrophil antibodies, 639
- Neutrophil antibody screen, 639
- Neutrophil gelatinase-associated lipocalin  
  (NGAL), 640–641, 641*t*
- Neutrophilia, 977–978*t*
- Neutrophils, 975  
  activated, lactoferrin from, 555

## 1048 index

- Newborn, metabolic screening for, 642–645
- Nexium, urea breath test and, 925
- Nicotinamide adenine dinucleotide (NADH) methemoglobin reductase enzyme, 623
- Nicotine, 646–648, 646*t*
- Nicotinic acid, catecholamine and, 956
- Nitrazine test strip, 51–52
- Nitrites, urinary, 935, 942
- Nitrofurantoin (Macrochantin, Nitrofan), urine color and, 932*t*
- Nitroglycerin  
catecholamines and, 956  
vanillylmandelic acid and, 955
- Nomogram, 9–10
- Non-A/non-B (NANB) hepatitis, 496
- Nonfunctioning nodule, in thyroid  
scanning, 883
- Nongranulocytes, 975
- Noninvasive prenatal paternity testing, 450
- Noninvasive prenatal testing (NIPT), 217–218
- Nonpancreatic diseases, amylase in, 56
- Nonpolyposis colorectal cancer, hereditary, 256
- Nonpregnant women, human chorionic gonadotropin in, 516
- Nonreactive fetal response, 415–416
- Nonsteroidal antiinflammatory drugs, thromboclastography and, 874–875
- Nonstrangulated ileocolic intussusception, barium enema in, 128
- Nonstress test (NST)  
fetal, 415–416  
fetal heart rate reactivity measured by, 406
- Norepinephrine, 954–957
- Normal serum K concentration, 724
- Normetanephrine, 954–957
- Normoblasts, in red blood cells, microscopic examination of, 150*b*
- Normochromic anemia, 774
- Normocytic cell, 772
- Nornicotine, 646–648
- Nose culture, 872–873  
NRAS, 588  
testing for, 588–589
- N-telopeptide (NTx), 168–169
- N-terminal fragment of pro-brain (B-type) natriuretic peptide (NT-pro-BNP), 636–637
- Nuclear imaging, salivary gland, 795–796
- Nuclear matrix protein 22 (NMP22), 145–146, 195–196*t*
- Nuclear medicine  
parathyroid scan and, 677  
salivary gland in, 795
- Nuclear scan, 1007  
carcinoid, 652–653  
cardiac, 207–209  
gallbladder, 431–432  
gallium, 433–434  
gastric emptying, 437–438, 437*t*  
gastrointestinal, 443–444, 444*t*  
for liver, 574  
of lung, 590  
Meckel diverticulum, 616–617  
neuroendocrine, 652–653  
octreotide scan, 652–653  
PET scan and, 719  
ProstaScint scan and, 739  
sentinel lymph node biopsy and, 805  
of spleen, 574  
thyroid, 883  
white blood cell, 979
- Nuclear stress test, 207–213  
myocardial, 207
- Nucleic acid amplification, for tuberculosis, 919–921. *See also* Tuberculosis (TB), testing
- Nucleosome (NCS), 68
- 5'-nucleotidase, 649, 649*t*
- Nystagmus, 357  
caloric study and, 192

## O

- Obesity, in abdominal ultrasound, 2
- Oblique view, in chest x-ray, 227
- Obstetric tests, 1013
- Obstetric ultrasonography, 685–687
- Obstruction series, 650, 651*t*  
of x-ray, 650
- Obstructive apnea, 827
- Obstructive biliary disease, 753
- Obstructive defects, pulmonary function tests and, 759
- Occlusal view, 322
- Occult blood, fecal, 844–846
- Occult disease, erythrocyte sedimentation rate for, 373
- Octreotide scan, 652–653
- Ocular photography, 424–425, 425*t*
- Oculovestibular reflex study, 192–193
- Odor, urinary, 931, 937, 940
- Oligoclonal gamma globulin bands, 578
- Oligohydrannios, 407
- Oligospermia, 801
- Oncogenes, mutations of, 588
- Oncology, PET scan and, 720–721
- Oncotype DX, 453–457  
in breast cancer, 454  
colon cancer assay, 256–257  
in prostate cancer, 744

Open lung biopsy, 586  
 Open pleural biopsy, 715  
 Ophthalmologic disease, fluorescein angiography for, 424  
 Opiates, 853*t*  
 Optical density, measurement of, in amniotic fluid, 50  
 Oral anticoagulant administration, 754  
 Oral contraceptives  
   activated partial thromboplastin time in, 682  
   testosterone and, 862  
   thyroxine-binding globulin and, 897  
   total thyroxine ( $T_4$ ) and, 895  
   triglycerides and, 908  
   triiodothyronine ( $T_3$  radioimmunoassay) and, 911  
 Oral glucose tolerance test (OGTT), 467–470, 468*f*  
 Oral mucosal transudate (OMT), in HIV-1 antibody tests, 508  
 Organ donors, anti-HTLV-I/II antibodies in, 525  
 Organ injury, 330  
 Organ-specific form, of autoimmune disease, 327  
 Orlistat, vitamin D and, 972  
 Oropharyngeal culture, 815  
 Oropharyngeal swabs, for SARS viral testing, 797  
 Osmolal gap, 654–655  
 Osmolality, 656–657  
   blood, 654–655, 655*t*  
   serum, 654–655, 655*t*  
   urine, 654, 656–657, 657*t*, 934  
 Osmosis, 371  
 Osmotic fragility (OF), 371–372, 372*t*  
 Osteocalcin, 168, 170  
 Osteoporosis  
   bone densitometry and, 157  
   bone turnover markers and, 169  
   vitamin D and, 971  
 O'Sullivan test, 465  
 Ova and parasites (O&P), stool for, 847–848  
 Ovarian cancer, 445–452, 446*t*  
   genetic testing in, 445–447  
 Overnight dexamethasone suppression test, 326  
 Oxidative hemolysis, Heinz bodies in, 481  
 Oximetry, 658–659, 659*f*, 659*t*  
   pulse, 658–659, 659*f*, 659*t*, 826  
 Oxycodone, 853*t*  
 Oxygen ( $O_2$ ) content, 106, 110  
 Oxygen ( $O_2$ ) saturation, 106, 109–110, 658–659, 659*f*, 659*t*  
 Oxygen ( $O_2$ ) saturation monitoring, 658  
   fetal, 417

Oxytocin challenge test (OCT), 409–411  
 Oxytocin stimulation, 415

## P

P wave, 342, 344*f*  
 p24 direct antigen, 506–509  
 p24 direct serologic antigen assay, 508  
 p53 protein, 176–178  
 Pacing, as stress testing, 211  
 Packed cell volume (PCV), 485–487  
 Packed red blood cell volume, 485–487  
 Pancreas, ultrasound of, 1–2, 4  
 Pancreatic carcinoma, amylase in, 56  
 Pancreatic duct flow, obstruction of, amylase in, 56  
 Pancreatic ducts, 364  
 Pancreatic elastase (PE), 660–661, 661*t*  
 Pancreatic enzymes, 662–664, 664*t*  
 Pancreatic hyperamylasemia, amylase in, 56  
 Pancreatic secretory test, 662–664, 664*t*  
 Pancreatic tests, 1013  
 Pancreatic tumor, gastrin-producing, 439  
 Pancreaticobiliary ducts, Manometric studies, 364  
 Pancreatitis  
   acute, lipase and, 562  
   amylase in, 56  
   chronic, 662  
 Pancreatobiliary FISH testing, 665–666  
 Pap smear, in HPV, 520–521  
 Pap test, 667–670, 669*f*, 670*t*  
 Papanicolaou (Pap) smear, 667–670, 669*f*, 670*t*  
 Papillary muscle, incision of, 364  
 PapSpin, 667–668  
 PapSure, 263–264  
 Paracentesis, 671–674, 672*f*, 673–674*t*  
 Paraneoplastic syndrome, antidiuretic hormone in, 72  
 Paraneoplastic tumors, prolactin and, 737  
 Parasites, in stool, 847  
 Parathormone, 675–676, 676*t*  
 Parathyroid hormone (PTH), 189–190, 675–676, 676*t*  
   vitamin D and, 971  
 Parathyroid scan, 677–678  
 Parathyroid scintigraphy, 677–678  
 Parathyroid tests, 1013  
 Parathyroids, 677  
 Parenchymal (hepatocellular) liver disease, 753–754  
 Parentage analysis, 445–452, 446*t*  
 Parietal cells, 90  
 Parkinson disease (PD), 679  
   PET scan and, 719–720  
   testing, 679–680

## 1050 index

- Parotid gland nuclear imaging, 795–796  
Paroxysmal nocturnal hemoglobinuria (PNH), 697  
PI-linked antigen, 697  
Partial pressure of O<sub>2</sub>, oximetry and, 658  
Partial thromboplastin time (PTT), 681–683, 683*t*  
Parvovirus B19, 684  
Parvovirus B19 antibody, 684, 684*t*  
Patau syndrome, 217  
Patch test, in allergy skin testing, 33, 35  
Paternity (parentage analysis), 445–452, 446*t*  
genetic testing of, 450  
PCO<sub>2</sub>, 106–107, 108*t*  
Peak inspiratory flow rate (PIFR), 761  
Peak level, drug, 332  
“Peak” test, in low molecular weight heparins, 77  
Pediatricians, IGF BPs in, 542  
Pelvic floor sphincter electromyography, 353–354  
for fecal incontinence, 353  
for urinary incontinence, 353  
Pelvic floor sphincter EMG, 353–354  
Pelvis  
computed tomography of, 272–277, 277*t*  
ultrasonography of, 685–687  
Penetrating peptic ulcer, amylase in, 56  
Penicillin  
in Bence Jones protein, 134  
cold agglutinins and, 253  
Pentagastrin stimulation, 187  
Pepsinogen, 688, 688*t*  
Pepsinogen II, 688  
Peptic ulceration, benign, 922  
Percutaneous IVUS imaging, 543  
Percutaneous lung biopsy, 843  
Percutaneous needle biopsy, 586, 715  
Percutaneous transhepatic cholangiography (PTHC), 364  
Percutaneous transluminal coronary angioplasty, 203  
Perfusion scan, 590  
cardiac, 207  
Periapical view, 322  
Pericardial fluid, pericardiocentesis and, 689, 690*f*  
Pericardiocentesis, 689–691, 690*f*, 691*t*  
Pericarditis, antimyocardial antibody in, 83  
Perinuclear ANCA (p-ANCA), 84–85  
Peripheral blood, PI-linked antigen, 697  
Peripheral blood smear, 149–151, 150*b*  
Peripheral synucleinopathy, 679  
Peripheral vascular atherosclerotic plaques, 770  
Peripheral vascular system, arteriography of, 116  
Periprosthetic joint infection (PJI), alpha defensin test in, 39  
Peritoneal fluid analysis, 671–674, 672*f*, 673–674*t*  
pericardiocentesis and, 671–674  
Peritoneal tap, 671–674, 672*f*, 673–674*t*  
Pernicious anemia (PA)  
intrinsic factor antibody in, 545  
pepsinogen and, 688  
vitamin B<sub>12</sub>, 969  
Persantine (dipyridamole), 211  
Pesticide exposure, cholinesterase and, 239  
PgR (progesterone receptor) assay, 735–736, 736*t*  
pH, 106–107, 108*t*  
esophageal, 380*f*  
fetal scalp blood, 417–418  
of pleural fluid, 867  
probe, acid reflux with, 377–378  
urinary, 931–933, 937–938, 941  
Phagocytosis, 975  
Pharmacogenetics, 334, 335*t*  
Phase-contrast magnetic resonance imaging (PC-MRI), 600  
Phenazopyridine (Pyridium)  
urinary bilirubin and, 939  
urine color and, 932*t*  
Phencyclidine, 853*t*  
Phenobarbital, 333*t*  
Phenolphthalein (Ex-Lax), urine color and, 932*t*  
Phenothiazines  
testosterone and, 862  
thromboelastography and, 874–875  
urine color and, 932*t*  
vanillylmandelic acid and, 955  
Phenotyping, alpha<sub>1</sub>-antitrypsin, 37–38, 38*t*  
Phenylketonuria, newborn metabolic screening and, 642–643  
Phenytoin, 333*t*  
free thyroxine (T<sub>4</sub>) and, 895  
thyroxine-binding globulin and, 897  
total thyroxine (T<sub>4</sub>) and, 895  
triiodothyronine (T<sub>3</sub> radioimmunoassay) and, 911  
urine color and, 932*t*  
vitamin D and, 972  
Pheochromocytoma, 692  
metanephrine, plasma free, 621  
suppression and provocative testing, 692–693  
vanillylmandelic acid and, 955  
phi (Prostate Health Index), 744  
Philadelphia chromosome, 195–196*t*  
Phlebography, 963–964  
Phosphate (PO<sub>4</sub>), 694–695, 695*t*  
Phosphatidylglycerol (PG), 50

- Phosphatidylinositol antigen, 697, 697*t*  
 Phosphatidylinositol glycan A (PIGA)  
   gene, 697  
 Phospholipase A2 receptor antibodies, 696  
 Phospholipid antibodies, 65  
 Phosphorus (P), 694–695, 695*t*  
   levels, 694  
 Photodiagnostic, oximetry, 658  
 Photography, ocular, 424–425, 425*t*  
 Photostimulation, 349  
 Physiologic jaundice, of newborn, 137  
*PIK3CA*, 588  
   testing for, 588–589  
 PI-linked antigen, 697, 697*t*  
 Pilocarpine iontophoresis, 857  
 Pinworms, 848  
 Pittsburgh Agent B (PIB), 58  
 Pituitary adenomas, prolactin level for, 737  
 Pituitary gland, anterior,  
   adrenocorticotrophic hormone in, 11  
 Pituitary reserve capacity, metyrapone test  
   in, 16  
 PLAC test, 564  
 Placenta, in human placental  
   lactogen, 523  
 Placental alpha microglobulin-1 (PAMG-1),  
   51–52  
 Placental grading, fetal biophysical profile  
   and, 407  
 Placental growth factor (PGF), 698, 698*t*  
 Plague, bioterrorism infectious agents  
   testing for, 142*t*, 143  
 Planes of reference, 343*f*  
 Plasma, porphyrin fractionation of, 717  
 Plasma ammonia levels, 47  
 Plasma anti-Xa assay, 77  
 Plasma coagulation system, 751  
 Plasma renin activity (PRA), 783–787, 787*t*  
 Plasma renin concentration (PRC),  
   783–787, 787*t*  
 Plasma testing, nicotine and metabolites  
   and, 647  
 Plasma thromboplastin antecedent, 252*t*  
 Plasmin, production of, 700  
 Plasminogen, 699, 699*t*  
   activators, 700  
 Plasminogen activator inhibitor 1 (PAI-1),  
   195–196*t*, 700–701, 701*t*  
   elevated levels of, 700  
 Plasminogen deficiency, 699  
 Platelet aggregation test, 702–703  
 Platelet agonist, 702  
 Platelet antibody detection, 704–705, 705*t*  
 Platelet closure time (PCT), 709–711, 710*t*  
 Platelet count, 706–708, 708*t*  
 Platelet function assay, 709–711, 711*t*  
 Platelet precursors, bone marrow biopsy in,  
   161–162  
 Platelets  
   dysfunction, 709  
   examination of, peripheral blood smear  
     in, 149  
   immune-mediated destruction  
     of, 704  
 Platelet-specific antigens, 704  
 Plethysmography, 713  
   arterial, 713–714  
 Pleural biopsy, 715–716  
 Pleural fluid analysis, 864–869, 869*t*  
   amylase in, 866  
   bacteriologic culture for, 866  
   carcinoembryonic antigen in, 867  
   cell counts in, 865  
   cultures for *Mycobacterium tuberculosis*  
     and fungus in, 866  
   cytology in, 867  
   explanation and related physiology in,  
     864–867  
   glucose in, 866  
   Gram stain for, 866  
   gross appearance of, 864–865  
   lactic dehydrogenase in, 866  
   pericardiocentesis and, 689  
   protein content in, 865–866  
   special tests for, 867  
   triglyceride in, 866  
 Pleural fluid tap, for SARS viral  
   testing, 797  
 Pleural tap, 864–869. *See also* Pleural fluid  
   analysis  
 Pneumoencephalography, brain computed  
   tomography and, 278  
 PO<sub>2</sub>, 106, 109  
 PO<sub>4</sub> (phosphate), 694–695, 695*t*  
 Poikilocytosis, in red blood cell, 149  
 Point-of-care home testing, 755  
 Poisoning, with carbon monoxide, 199  
 Polyclonal gamma globulin  
   elevation, 748*t*  
 Polycystic ovary syndrome, adrenal steroid  
   precursors in, 19  
 Polycythemia, 375  
 Polycythemia vera, 770  
 Polyhydramnios, 407–408  
 Polymerase chain reaction (PCR) culture,  
   for tuberculosis, 917  
 Polymorphonuclear leukocytes (PMNs),  
   975  
 Polymorphonucleocyte antibodies (PMN  
   ab), 639  
 Polys (polymorphonuclear leukocytes),  
   975  
 Polysomic cells, 665  
 Polysomnography (PSG), 826–829  
 Pontiac fever, 560  
 Porphobilinogen deaminase, 953

## 1052 index

- Porphobilinogens, 717–718, 717–718*t*  
Porphyria, 321  
  porphyrins and porphobilinogens, 717  
  uroporphyrinogen-1-synthase and, 953  
Porphyrins, 717–718, 717–718*t*  
Portal hypertension, 574, 770  
Positron emission mammography (PEM), 721  
Positron emission tomography (PET), 433, 719–723, 723*t*  
  with amyloid, for Alzheimer disease, 58  
Postcoital cervical mucus test, 822–823  
Postcoital test, 822–823  
Posteroanterior (PA) view, in chest x-ray, 227  
Postinfarction mural thrombi, echocardiography for, 339  
Postmortem tissue, for SARS viral testing, 797  
Postmyocardial infarction, antimyocardial antibody in, 83  
Postprandial glucose (PPG), 465–466, 466*t*  
Posttransfusion purpura, 704  
Potassium  
  blood, 724–726  
  with calcium, and magnesium, 597  
  urine, 727, 727*t*  
Potassium iodide, thyroid-stimulating hormone and, 887  
PR interval, 342, 344*f*  
Prealbumin (PAB), 728–729, 729*t*  
Preeclampsia, 698  
Pregnancy  
  activated partial thromboplastin time in, 682  
  alpha-fetoprotein in, 41  
  amino acids in, 45  
  antithrombin III in, 97  
  ectopic, human chorionic gonadotropin in, 515–516  
  folic acid in, 426  
  high-risk, pregnanediol, 731  
  human placental lactogen in, 523  
  morbidity, anticardiolipin antibodies in, 65  
  pelvic ultrasonography in, 685–687  
  plasminogen and, 699  
  progesterone in, 733  
  rubella antibody test in, 793–794  
  toxoplasmosis antibody titer for, 902  
  transferrin in, 548  
Pregnancy tests  
  human chorionic gonadotropin in, 515–517  
  sexual assault testing and, 810  
Pregnancy-associated plasma protein-A (PAPP-A), 730  
Pregnanediol, 731–732, 732*t*  
Pregnenolone, 18–19, 18–19*t*  
Pre-mature rupture, of membranes, 51–52  
Premedication, before EMG, 351  
Prenatal cell-free DNA screening, 217–218  
Prenatal testing, noninvasive, 217–218  
Prenatal tests, 1013  
Preoperative parathyroid scanning, 677  
Prerenal azotemia, 837  
Pressure, of cerebrospinal fluid, 577  
Presurgery, total blood volume for, 900  
Prevacid, urea breath test and, 925  
Prick-puncture test, in allergy skin testing, 33, 35  
Prilosec, urea breath test and, 925  
Primary adrenal gland failure, adrenocorticotropic hormone in, 11  
Primary adrenal insufficiency  
  adrenocorticotropic hormone stimulation test with cosyntropin in, 13–14  
  chronic, 21-hydroxylase antibodies in, 530  
Primary aldosteronism, aldosterone in, 25–28  
Primary biliary cirrhosis, 81  
  antimitochondrial antibody in, 82  
Primary constriction, centromere in, 67  
Primary hypothyroidism, 887*t*, 888  
Primary liver diseases, haptoglobin in, 479  
Primary membranous nephropathy (pMN), 696  
Primary muscular disorders, aldolase in, 23  
Primary renal diseases, antidiuretic hormone in, 71–72  
Primary skin disorders, fetal skin biopsies for, 419  
Proaccelerin, 249–250*t*, 252*t*  
Proconvertin, 249–250*t*, 252*t*  
Proctoscopy, 258–261  
Progesterone  
  endometrial biopsy and, 362  
  major effect of, 733  
  receptor assay, 735–736, 736*t*  
Progesterone assay, 733–734, 734*t*  
Progesterone receptor (PR) assay, 735–736, 736*t*  
Programmed death ligand 1 (PD-L1), 195–196*t*  
Progressive systemic sclerosis (PSS), antiscleroderma antibody in, 91  
Proinsulin C-peptide, 293–294, 294*t*  
Prolactin levels (PRLs), 737–738, 738*t*  
Prolactin stimulation tests, 737  
Prolactin suppression tests, 737  
Prolactin-secreting pituitary acidophilic adenoma, 737  
Proliferative-type endometrium, 362  
Prolonged test, dexamethasone suppression test, 324–326  
Promark testing, in prostate cancer, 744  
Propoxyphene, 853*t*



- Propranolol, 333*t*  
   free thyroxine ( $T_4$ ) and, 895  
   thyroxine-binding globulin and, 897  
   total thyroxine ( $T_4$ ) and, 895  
   triiodothyronine ( $T_3$  radioimmunoassay)  
   and, 911
- Propylthiouracil, total thyroxine ( $T_4$ ) and,  
 895
- Prostaglandin D(2) (PGD[2]), 136
- ProstaScint scan, 739–740
- Prostate cancer, 739  
   acid phosphatase in, 7  
   genomic testing of, 454–455  
   metastatic, 739  
   PSA in, 743–745, 745*t*
- Prostate cancer gene 3 (PCA3), 743–745,  
 745*t*
- Prostate cancer-specific biomarkers, 744
- Prostate disease, PSA and, 743
- Prostate gland  
   disease outside, 739  
   ultrasound of, 741–742
- Prostate Health Index (phi), 744
- Prostate/rectal sonogram, 741–742, 741*f*
- Prostate-specific antigen (PSA), 195–196*t*,  
 741–745, 745*t*  
   age-adjusted, 743  
   density, 743  
   free, 743  
   isoforms, 744  
   velocity, 743
- Prostate-specific proteins, 743
- Prostatic acid phosphatase (PAP), 7–8, 8*t*
- Prostatic-specific membrane antigen, 743
- Prostatitis, PSA and, 743
- Protamine, nomogram in, 9–10
- Protamine sulfate, 78
- Protein, 746–750, 750*t*  
   activated, 397  
   acute-phase, 373  
   acute-phase reactant, 404  
   ammonia and, 47  
   anion gap and, 63  
   blood, 746  
   in cerebrospinal fluid, 578  
   peritoneal fluid and, 671–672  
   PGF and, 698  
   in pleural fluid analysis, 865–866  
   PSA and, 743  
   total, 746–750, 750*t*  
   urinary, 746, 933, 938, 941
- Protein C, 751–752, 752*t*  
   activated, 682  
   endogenous anticoagulant, 397
- Protein electrophoresis, 746–750, 748*t*, 750*t*  
   for AAT deficiency, 37  
   of urine, 134
- Protein S, 751–752, 752*t*
- Proteinase 3 (PR3) autoantigen, 84
- Proteinase 3-ANCA (PR3-ANCA), 84
- Protein/creatinine ratio, 933
- Proteinuria, 933
- Proteus vulgaris* antigens, 399
- Prothrombin, 248, 249–250*t*, 252*t*
- Prothrombin fragment (F1 + 2), 877, 878*t*
- Prothrombin time (PT), 753–756, 756*t*
- Protidine, 753–756, 756*t*
- Proton pump inhibitors, urea breath test  
   and, 925
- Protonix, urea breath test and, 925
- Protoporphyrin, 987
- Provocative tests, pheochromocytoma  
   suppression and provocative testing,  
   692
- Pseudocholine challenge, 763
- Pseudocholinesterase (PChE), 238–240,  
 240*t*
- Pseudohyperparathyroidism, 675
- Pseudomonas* organisms, 931
- Pulmonary angiography, 757–758
- Pulmonary anthrax, 141
- Pulmonary arteriography, 757–758
- Pulmonary artery pressure, 202–203*t*
- Pulmonary embolism (PE)  
   D-dimer, 319  
   lung scan, 590  
   pulmonary arteriography and, 757
- Pulmonary embolus, ischemia-modified  
   albumin in, 550
- Pulmonary function tests (PFTs), 759–763,  
 763*t*
- Pulmonary scintiphotography, 590–592  
   contraindications to, 591
- Pulmonary system, tests for, 996–997
- Pulmonary wedge pressure, 202–203*t*
- Pulse oximetry, 110, 658–659, 659*f*, 659*t*,  
 826
- Punch biopsy, 224
- Purified protein derivative (PPD) test,  
 915–916
- Purpura  
   cryoglobulin and, 306  
   idiopathic thrombocytopenia, 704  
   posttransfusion, 704
- Pyclography, 764–768
- Pyocelle, 799
- Pyridinium (PYD) crosslinks, 170
- Pyridium, urine, 168
- Q**
- QRS complex, 342, 344*f*
- Quadruple test, 612
- Qualitative fetal hemoglobin stool test, 105
- QuantiFERON-TB Gold (QFT-G test),  
 919–921. *See also* Tuberculosis (TB),  
 testing

## 1054 index

- Quantitative computed tomography (QCT),  
in bone mineral density, 157
- Quantitative fibrinogen, 422–423, 423 $t$
- Quantitative stool fat determination,  
402–403, 403 $t$
- Quantitative ultrasound, in bone mineral  
density, 157
- Queckenstedt-Stookey test, 581–582
- Quinidine, catecholamine and, 956
- Quinine, WBC levels and, 976
- R**
- R factor, 63–64, 64 $t$
- Rabies, 769
- Rabies-neutralizing antibody test, 769
- Radioactive carbon ( $^{13}\text{C}$ ), in breath test, 483
- Radioactive chemicals, 719
- Radioallergosorbent test (RAST), 31–32,  
32 $t$
- Radioimmunoscintigraphy (RIS), 739–740
- Radioiodine imaging, in thyroglobulin, 880
- Radioiodine-labeled albumin, 900
- Radionuclear scan, 780
- Radionuclide, in bone scan, 165
- Radionuclide renal imaging, 779–782
- Radiopaque dye, in subarachnoid space, 632
- Radiorenography, 779–782
- Rape testing, 809–812
- Rapid plasma reagin (RPR), 859–860
- Rapid stimulation test, 13
- Rapid test  
in adrenocorticotrophic hormone  
stimulation test with cosyntropin, 14  
dexamethasone suppression test, 326  
mononucleosis, 630
- Rapid urease testing, for *Helicobacter  
pylori*, 482
- Raynaud phenomenon, cryoglobulin and,  
306
- Reactant protein, 373
- Reactive fetal response, 415
- Reagin, 859
- Recombinant human erythropoietins,  
transferrin receptor (TfR) assay  
and, 906
- Rectal culture, 816 $f$
- Rectal EMG procedure, 353–354  
for fecal incontinence, 353  
for urinary incontinence, 353
- Rectal ultrasonography, 741–742, 741 $f$
- Red blood cells (RBCs), 371  
in cerebrospinal fluid, 578  
fetal, leakage of, 413  
fragility, 371–372  
group of, 152  
Heinz bodies in, 481  
peritoneal fluid and, 671–672
- Red blood cells (RBCs) (*Continued*)  
in pleural fluid analysis, 865  
sedimentation rate, 373–374  
smear, 149–151, 150 $b$   
urinary, 937, 939, 944  
values, 770  
zinc protoporphyrin in, 987
- Red blood cell (RBC) casts, 937
- Red blood cell (RBC) cholinesterase,  
238–240, 240 $t$
- Red blood cell (RBC) count, 770–771,  
771 $t$
- Red blood cell (RBC) distribution width  
(RDW), 774
- Red blood cell (RBC) indices, 772–775,  
773 $t$   
abnormal findings in, 775, 775 $t$   
interfering factors in, 774
- Red blood cell (RBC) volume, 900–901
- Reductase, 935
- Reflectance meter, 464
- Reflux, of barium, 128, 131
- Regadenoson, 207
- Renal angiography, 114
- Renal biopsy, 776–778, 777 $f$   
in Goodpasture syndrome, 79
- Renal blood flow (perfusion) scan, 779
- Renal calculus analysis, 945–946
- Renal failure  
in aluminum, 43  
creatinine and, 301  
creatinine clearance and, 304  
lipase and, 562  
neutrophil gelatinase-associated lipocalin  
for, 640
- Renal function scans, 780
- Renal function studies, 301  
blood urea nitrogen in, 155
- Renal hypertension scan, 780
- Renal obstruction scan, 780
- Renal panel, 1013
- Renal scanning, 779–782
- Renal structural scan, 779
- Renal system, tests for, 997–998
- Renal vein, 783–787  
assays, 784
- Renin, 60
- Renin assay, 783–787, 787 $t$
- Renin stimulation test, 784–785
- Renin-angiotensin system, 783
- Renogram, 780
- Renovascular hypertension, 783, 784 $f$
- Reproductive system, tests for, 998–999
- Reserpine  
catecholamines and, 956  
triiodothyronine ( $\text{T}_3$  radioimmunoassay)  
and, 911  
vanillylmandelic acid and, 955

Residual volume (RV), 760*f*, 761  
 Resin, in blood culture, 147  
 Respiratory acidosis, 108–109*t*  
 Respiratory alkalosis, 108–109*t*  
 Respiratory distress syndrome (RDS), in  
 L/S ratio, 49–50  
 Respiratory disturbance index (RDI), in  
 sleep studies, 826  
 Respiratory specimens, types of, for SARS  
 viral testing, 797  
 Restrictive defects, pulmonary function tests  
 and, 759  
 RET, 588  
 testing for, 588–589  
 RET germline mutations, 449  
 Retic count, 788–789, 789*t*  
 Reticulocyte, 788  
 Reticulocyte count, 788–789, 789*t*  
 Reticulocyte hemoglobin equivalent,  
 788–789  
 Reticulocyte index, 788  
 Reticulocyte-specific hemoglobin content,  
 788–789  
 Retinol binding protein, 628–629  
 Retrograde cholangiopancreatography,  
 endoscopic, 364–367, 366*f*, 665  
 for jaundice, 364  
 Retrograde pycnography, 764–768  
 Rh factor, in blood typing, 153  
 Rh isoimmunization, fetal status affected  
 by, 51  
 Rheumatic diseases, antinuclear antibody  
 in, 86  
 Rheumatoid arthritis (RA), 790  
 anticyclic citrullinated peptide antibody  
 in, 70  
 Rheumatoid factor (RF), 790–791, 791*t*  
 RhoGAM, 153, 421  
 Riboflavin  
 catecholamines and, 956  
 urine color and, 932*t*  
 Ribosomal P Ab, 792, 792*t*  
 Ribosome, 792  
 Ribosome P antibodies, 792, 792*t*  
 Rickettsial species, 399  
 Rifampicin, free thyroxine (T<sub>4</sub>) and, 895  
 Rifampin, urine color and, 932*t*  
 Right-sided heart catheterization, 203  
 Ristocetin, 702  
 test, 702–703  
 Rivaroxaban, 77  
 RNA polymerase III antibody, 91  
 RNA quantification, of HIV, 503–505  
 RNA testing, HCV, 496  
 Romberg sign, 192  
 ROME III criteria, 316  
 ROSI, 588  
 testing for, 588–589

Round cells (spherocytes), in osmotic  
 fragility, 371  
 Routine blood pressure, 202–203*t*  
 Routine serum protein electrophoresis, for  
 AAT deficiency, 37  
 Rubber stoppers, aluminum silicate in, 43  
 Rubella, 793  
 Rubella antibody test, 793–794, 793*t*

**S**

Sacral x-ray studies, 839–840  
 Salicylates, 333*t*  
 catecholamines and, 956  
 triiodothyronine (T<sub>3</sub> radioimmunoassay)  
 and, 911  
 Saline infusion sonography (SIS), 686  
 Saline-extracted antigen, 75  
 Saliva  
 cortisol in, 290–292, 292*t*  
 for drug testing, 852  
 Salivary cortisol, 290–292, 292*t*  
 Salivary duct, x-ray of, 818  
 Salivary gland nuclear imaging, 795–796  
 Salivary hyperamylasemia, amylase in, 56  
 Sandhoff disease, 500  
 Sarcoidosis, angiotensin-converting enzyme  
 in, 61  
 Scalp blood pH, fetal, 417–418  
 Scintigraphy. *See also* Nuclear scan  
 abdominal, 443–444, 444*t*  
 gastrointestinal, 443–444, 444*t*  
 hepatobiliary, 431–432  
 Scintillator detector/camera, in bone  
 densitometry, 159  
 Scintiphography, pulmonary, 590–592  
 contraindications to, 591  
 Scl-70 antibody, 91  
 Scleroderma, 91, 790  
 Scleroderma antibody, 91  
 Scratch test, in allergy skin testing, 33  
 Scrotal ultrasound, 799–800, 800*t*  
 Secondary adrenal insufficiency  
 adrenocorticotrophic hormone stimulation  
 test with cosyntropin in, 13–14  
 metyrapone test in, 16  
 Secondary aldosteronism, aldosterone in,  
 25, 27–28  
 Secondary hypothyroidism, 886, 887*t*, 888  
 Second-strength PPD, 915  
 Secretin and pancreozymin, 662  
 Secretin test, 440  
 Secretory-type endometrium, 362  
 Sed rate test, 373–374, 374*t*  
 Sedation, before EMG, 351  
 Seizure disorder, electroencephalography  
 for, 349  
 Semen, antisperm antibodies in, 93

## 1056 index

- Semen analysis, 801–804
- Seminal fluid, acid phosphatase in, 7
- Sentinel lymph node biopsy (SLNB), 805–806
- Sequential multiple analysis (SMA) 7, 133
- Sequential multiple analysis (SMA) 12, 270–271, 270*t*
- Sequential multiple analysis (SMA) 20, 270–271, 270*t*
- Serine protease inhibitor, 97–98, 98*t*
- Seroconversion window, 506
- Serologic test for syphilis (STS), 859–860
- Serologic tests
- for *Helicobacter pylori*, 483
  - for herpes simplex, 498
  - blood for, 499
- Serotonin, 807–808, 808*t*
- Serum
- antispermat antibodies in, 93
  - nicotine and metabolites for, 646*t*, 647
  - for SARS viral testing, 797
- Serum albumin, 747
- Serum aldolase, 23
- Serum amylase, 56
- Serum angiotensin-converting enzyme (SACE), 61–62, 61–62*t*
- Serum calcium, 189–191, 189*t*, 191*t*
- Serum chromogranin A, serotonin and, 807
- Serum CO<sub>2</sub> test, 197
- Serum cortisol, 290–292, 292*t*
- Serum creatinine, 301–302, 302*t*
- Serum enzyme tests, before EMG, 350
- Serum ferritin study, 404
- Serum glutamic-oxaloacetic transaminase (SGOT), 125–127, 126–127*t*
- Serum glutamic-pyruvic transaminase (SGPT), 21–22, 22*t*
- Serum haptoglobin test, 479
- Serum hepatitis, 493–494
- Serum iron, 546–547
- Serum osmolality, 654–655, 655*t*
- test, 654
- Serum progesterone level, progesterone assay and, 733
- Serum protein electrophoresis (SPEP), 746–750, 750*t*, 754*t*
- for AAT deficiency, 37
- Serum protein quantification, 536
- Sestamibi cardiac scan, 207–209
- Severe acute respiratory syndrome (SARS) viral testing, 797–798
- Sex, of fetus, 51
- Sexual assault evidence collection kit, 809
- Sexual assault testing, 809–812
- Sexually transmitted disease testing, 813–817
- Sexually transmitted diseases (STDs), 813*t*
- Chlamydia* infection as, 230
  - Shake test, 50
  - Shift to the left, in WBC production, 975
  - Short-gut syndrome, 402
  - Shortness of breath, 636
  - Sialography, 818–819
  - Sickle cell anemia, 820
    - newborn metabolic screening and, 643  - Sickle cell disease, sickle cell screen in, 820
  - Sickle cell screen, 820–821
  - Sickle cell trait, 820
  - Sickledex, 820–821
  - Sigmoidoscopy, 258–261
  - Signal-averaged ECG (SAECG), 344
  - Signaling proteins, activation of, 588
  - Simple cervical biopsy, 224
  - Sims-Huhner test, 822–823
  - Single contrast barium enema, 128
  - Single tracer double phase (STDTP) method, 677–678
  - Single-photon emission computed tomography (SPECT), 433, 679
    - cardiac, 207
    - of liver, 574, 581*f*  - Single-stranded DNA, antibody against, 74
  - Sjögren antibodies, 95–96
  - Sjögren syndrome, 95, 790
  - Skeletal muscle, disease or trauma of, and myoglobin, 635
  - Skeletal system, tests for, 999
  - Skin biopsy, 824
  - Skin biopsy antibodies, 824
  - Skin immunohistopathology, 824
  - Skin testing, allergy, 33–36, 36*t*
    - in IgE antibodies, 31
    - intradermal method in, 35
    - patch method in, 35
    - prick-puncture method in, 35  - Skinny-needle thyroid biopsy, 881–882, 881*t*
  - Skull x-ray, 825
  - Sleep apnea, sleep studies and, 827
  - Sleep disorders, 826
  - Sleep EEG, 349
  - Sleep studies, 826–829
  - SMA 20, 270–271, 270*t*
  - SMAC 7, 133
  - Small bowel enema, 830–832
  - Small bowel follow-through (SBF), 830–832
  - Small cell lung cancer (SCLC), 638
  - Small intestinal bacterial overgrowth (SIBO) test, 833–834
  - Smallpox, bioterrorism infectious agents testing for, 142*t*, 143
  - Smears, sexually transmitted disease, 815
  - Society for Maternal-Fetal Medicine (SMFM), 217
  - Sodium fluoride F18 injection (18F NaF) scan, PET scan and, 721

- Sodium (Na)  
 in blood, 835–836, 836*t*  
 fractional excretion of, 837  
 in urine, 837–838, 838*t*  
 values in children, 857
- Sodium resorption, 724
- Soft tissue swelling, in bone x-ray, 172
- Soluble fms-like tyrosine kinase-1 (sFlt-1), 698
- Somatomedin C, 541–542, 541–542*t*
- Somatomedins, 474, 541
- Somatosensory-evoked responses (SERs), 394  
 abnormal latency for, 396
- Somatotropin hormone (SH), 474–476, 476*t*
- Sonogram, breast, 181–182, 182*t*
- Sound sensors, 826
- Specific gravity, urinary, 934–935, 938, 942
- Specimen collection, Pap test, 667–668
- Spectroscopy, 263–264
- Speculoscopy, 263–264
- Sperm, agglutination and inhibition of, 93–94
- Sperm chromatin structure assay test, 802
- Sperm count, 801–804
- Sperm DNA fragmentation assay (SDFA) test, 802
- Sperm DNA integrity, 802
- Sperm examination, 801–804
- Sperm penetration assay (SPA), 802
- Sperma antibodies, 93–94
- Sperm-cervical mucus interaction, 802
- S-phase fraction, 176–178
- Spherocytes (round cells)  
 in osmotic fragility, 371  
 in red blood cells, microscopic examination of, 150*b*
- Spherocytosis  
 Heinz bodies in, 481  
 hereditary, erythrocyte fragility for, 371
- Sphincter of Oddi, Manometric studies of, 364
- Spinal tap, 576–583. *See also* Cerebrospinal fluid
- Spinal x-ray, 839–840
- Spine, cervical or lumbar, MRI of, 601
- Spinnbarkeit, 822
- Spirometry, 759, 761–762
- Spironolactone, 785  
 testosterone and, 862
- Spleen  
 enlarged, 770  
 scanning of, 574–575  
 contraindications to, 574
- Split study, sleep and, 827
- Spontaneous muscle movement, during EMG, 350
- Sputum  
 collection, in tuberculosis culture, 918  
 for SARS viral testing, 797
- Sputum culture and sensitivity, 841–842
- Sputum cytology, 843
- Sputum Gram stain, 841–842
- Sputum smear, 917
- Sputum tests, 1008
- Squamocolumnar junction, 669
- Squamous cell carcinoma (SCC) antigen, 195–196*t*
- ST segment, 343, 344*f*
- Stab cells, 975
- Stable factor, 249–250*t*, 252*t*
- Staining, of sputum, 841
- Standard deviation (SD), in bone mineral density, 157–158
- Standardized uptake value (SUV), 720–721
- Statins, triglycerides and, 908
- STD (sexually transmitted disease) culture, 814
- STD testing, 813–817
- Steatorrhea, 402
- Stein- Leventhal syndrome, adrenal steroid precursors in, 19
- Stenosis, of coronary artery, intravascular ultrasound in, 543
- Stent restenosis, intravascular ultrasound in, 543
- Stenting, 114
- Stents  
 ERCP and, 364  
 intracoronary, 203
- Stereotactic biopsy, 606, 607*f*
- Steroids  
 anabolic  
 thyroxine-binding globulin and, 897  
 total thyroxine ( $T_4$ ) and, 895  
 testosterone and, 862  
 thyroid-stimulating hormone and, 887  
 thyroxine-binding globulin and, 897  
 WBC levels and, 976
- Stool  
 culture and sensitivity (C&S) of, 847–848  
 DNA home test of, 846  
 DNA sample of, 844–846  
 fat in, 402, 403*t*  
 Hemocult slide test of, 846  
 for leukocytes, 849  
 for occult blood (OB) testing, 844  
 for ova and parasites, 847–848  
 pancreatic elastase in, 660–661  
 for swallowed blood, 105  
 tablet test of, 846
- Stool cancer screening, 844–846
- Stool culture, 847–848
- Stool fat determination, quantitative, 402–403, 403*t*

## 1058 index

- Stool tests, 1008
- Strep screens, throat and nose cultures for, 872
- Streptococcal infection, throat and nose cultures for, 872
- Streptococcus* group B antigen detection, 850–851, 851*t*
- Streptococcus* serologic testing, 850–851, 851*t*
- Streptokinase, thrombosis indicators fibrin monomers and, 877
- Streptolysin O, 850–851
- Streptozyme, 850–851, 851*t*
- Stress testing, 210
- cardiac, 210–213, 339–340
  - chemical, 210
  - echo, 210–213
  - exercise, 210–213
- Striated muscle antibody, IgG, 5
- Stroke, 719–720
- acute, ischemia-modified albumin in, 550
- Stroke volume (SV), 202–203*t*
- Strongyloides* spp. (tapeworm), in stool, 847
- Structural renal scans, 781
- Stuart factor, 249–250*t*, 252*t*
- Subarachnoid space, pressure in, measurement of, 577
- Substance abuse testing, 852–854, 853*t*
- Subxiphoid route, pericardiocentesis using, 690*f*
- Succinylcholine, 238
- Sucralfate (Carafate), urea breath test and, 925
- Sulfasalazine (Azulfidine), urine color and, 932*t*
- Sulfonamides, WBC levels and, 976
- Superconducting quantum interference device (SQID), 347
- Superficial bladder cancers, recurrence rate for, 145
- Superscan, in bone scan, 166
- Supine abdominal x-ray study, 650
- Suppression test
- ACTH, 324–326, 326*t*
  - cortisol, 324–326, 326*t*
  - dexamethasone, 324–326
  - overnight, 326
- Suprapubic aspiration, of urine, 948
- Surfactant activity, measurement of, 50
- Surgical tubal ligation, hysterosalpingography in, 531
- Swallowed blood, stool for, 105
- Swallowing examination, 855–856
- Sweat, for drug testing, 852
- Sweat electrolytes test, 857–858
- Symbols and units of measurement, 1015–1016, 1015–1016*t*
- Symptomatic ulcer disease, *Helicobacter pylori* testing in, 482
- Syndrome of inappropriate antidiuretic hormone (SIADH), antidiuretic hormone in, 72
- Synovasure, 39–40
- Synovial fluid, 118
- analysis, arthrocentesis with, 118–121, 121*t*
- Synovial fluid glucose, 119
- Synucleinopathy, 679
- Syphilis, 813*t*
- serology for, in cerebrospinal fluid, 579
  - in sexually transmitted disease testing, 813–814
- Syphilis detection test, 859–860
- Systemic lupus erythematosus (SLE), 790
- anti-DNA test in, 74
  - antiextractable nuclear antigens in, 75
  - antinuclear antibody (ANA) for, 86, 86*t*
  - antinucleosome antibody in, 68
- Systemic mast cell disease (SMCD), 11 beta-prostaglandin F(2) alpha in, 136
- Systemic mastocytosis, 11 beta-prostaglandin F(2) alpha in, 136
- Systolic left ventricular pressure, 202–203*t*

## T

- T wave, 343, 344*f*
- T<sub>3</sub>, thyroid-stimulating hormone and, 887
- T<sub>3</sub> radioimmunoassay, 910–911, 911*t*
- T<sub>3</sub> toxicosis, 910
- Tablet test, of stool, 846
- Tamoxifen, 686
- thyroxine-binding globulin and, 897
- Tandem mass spectrometry, 642
- Tap test, 50
- Tape test, 848
- Target heart rate, 210–211
- Tartrate-resistant acid phosphatase (TRAP), 7–8, 8*t*
- Tau protein, 58
- Tay-Sachs disease (TSD), 445–452, 446*t*, 500
- genetic testing of, 447–448
- Tc 99m<sup>99m</sup>, in thyroid scanning, 883
- T-cells, 975
- Technetium agents, 207–208
- Technetium-99 m (<sup>99m</sup>Tc), 431
- Technetium-99 m diethylenetriamine pentaacetic acid (<sup>99m</sup>Tc-DTPA), 780
- Technetium-99m-mercaptoacetyltryglycine (<sup>99m</sup>Tc-MAG3), 780
- Temperature regulation, febrile antibodies and, 399
- Terminal ileum, barium in, 128

- “Test meal”, 437
- Testes, ultrasound of, 799–800, 800*t*
- Testosterone, 861, 861*t*, 863*t*
  - free, 861, 861*t*
  - total, 861*t*
- Testosterone stimulation tests, 862
- Tetracyclines, catecholamines and, 956
- Thalassemia, erythrocyte fragility for, 371
- Thallium scan, 207–209
- T-helper cells, 219
- Theophylline, 333*t*
  - thromboelastography and, 874–875
- Therapeutic drug monitoring (TDM), 332–335, 333*t*
- ThinPrep Imaging System, for Pap smear, 667–670, 670*t*
- Thiopurine methyltransferase (TPMT), 334
- Thoracentesis, 864–869, 865*f*, 869*t*. *See also*
  - Pleural fluid analysis
  - pleural biopsy and, 715
- Thoracic x-ray studies, 839–840
- Thoroscopic lung biopsy, 586
- Thoracoscopy, 870–871
- Thoracotomy, video-assisted (VAT), 870
- 3-day test, in adrenocorticotrophic hormone stimulation test with cosyntropin, 14
- Three-dimensional echocardiography, 339
- Three-dimensional mammography, 606–607
- Three-phase bone scan, 165–166
- Throat culture, 872–873
- Thrombocyte count, 706–708, 708*t*
- Thrombocytopenia, 706
- Thrombocytopenia, 704, 706, 712
  - causes of, 706
- Thrombocytopenic disorders, differential diagnosis of, 712
- Thrombocytosis, 706
- Thromboelastography, 874–876
- Thromboelastometry, 874–876
- Thrombophilia, FVL-associated factor V Leiden for, 397
- Thromboplastin, 252*t*
- Thromboplastin time, partial/activated, 671–674, 672*f*, 673–674*t*
- Thrombosis indicators fibrin monomers, 877–878, 878*t*
- Thromboxane A<sub>2</sub>, 710
- Thyretin, 728–729, 729*t*
- Thyrocalcitonin, 187–188, 188*t*
- Thyrogen, 879–880
- Thyrogen-stimulated testing, 879
- Thyrogen-stimulated thyroglobulin, 879–880, 880*t*
- Thyroglobulin antibody, 99, 99*t*
- Thyroglobulin (Tg), 195–196*t*, 879–880, 879–880*t*
- Thyroid antithyroglobulin antibody, 99, 99*t*
- Thyroid autoantibody, 99–100, 99–100*t*
- Thyroid cancer, 445–452, 446*t*, 879
  - genetic testing of, 449–450, 455–456
- Thyroid diseases, antithyroid peroxidase antibody in, 100
- Thyroid echogram, 891
- Thyroid fine needle aspiration biopsy (FNAB), 881–882, 881*t*
- Thyroid lobes, 677
- Thyroid nodules, 883
- Thyroid scanning, 883–885
- Thyroid scintiscan, 883–885
- Thyroid screening panel, 1014
- Thyroid sonogram, 891
- Thyroid ultrasound, 891, 891*t*
- Thyroid-binding inhibitory immunoglobulin (TBII), 889–890
- Thyroiditis
  - antithyroglobulin antibody in, 99
  - antithyroid peroxidase antibody in, 100
- Thyroid-releasing hormone (TRH) stimulation test, 886–887
- Thyroid-stimulating hormone receptor (TSHR) antibodies, 889–890
- Thyroid-stimulating hormone (TSH), 879, 886–887, 887*t*
- Thyroid-stimulating hormone (TSH) stimulation test, 888
- Thyroid-stimulating immunoglobulins (TSIs), 889–890, 890*t*
- Thyrotropin, 886–887
- Thyrotropin receptor antibody, 889–890
- Thyrotropin-releasing factor (TRF) stimulation test, 892–893
- Thyrotropin-releasing hormone (TRH) stimulation test, 892–893
- Thyroxine (T<sub>4</sub>), 894–896
  - abnormal findings in, 896, 896*t*
  - interfering factors in, 895
  - screen, 894–896
- Thyroxine-binding globulin (TBG), 897–898, 897–898*t*
- Thyroxine-binding prealbumin (TBPA), 728–729, 729*t*
- Tidal volume (TV), 760*f*, 761
- Tilt-table test, 359, 361
- Tissue, in antiglomerular basement membrane antibodies, 79
- Tissue factor, 252*t*
- Tissue folate, 426
- Tissue plasminogen activator (tPA), 700
- Tissue transglutaminase antibodies (tTG-Ab), 458–459, 458*t*
- Tobramycin, 333*t*
- Toluidine blue dye test, 810
- TORCH antibody panel, 1014

## 1060 index

- TORCH (toxoplasmosis, other, rubella, cytomegalovirus, herpes) test, 793–794, 899  
  cytomegalovirus, 317
- Total blood volume (TBV), 900–901, 901*t*
- Total CD4 count, 219
- Total complement, 266–268, 267–268*t*
- Total cortisol, 290–292, 292*t*
- Total hexosaminidase, 500–501
- Total IgE serum, 31
- Total iron-binding capacity (TIBC), 546–549, 549*t*  
  transferrin and, 547  
  transferrin saturation and, 547–548
- Total lung capacity (TLC), 760*f*, 761
- Total parenteral nutrition (TPN), 728
- Total serum bilirubin level, 138–139
- Total thyroxine, 894–896. *See also* Thyroxine (T<sub>4</sub>)
- Total/ionized calcium, 189–191, 189*f*, 191*t*
- Toxicology screening, 852–854, 853*t*  
  for urine, 1014
- Toxin B, cytolethal distending, 316
- Toxin B gene, *C. difficile*, 246
- Toxin C gene, *C. difficile*, 246
- Toxoplasma gondii*, 902
- Toxoplasmosis, 902  
  newborn metabolic screening and, 644  
  TORCH test for, 899
- Toxoplasmosis antibody titer, 902
- Tracheal aspirate, for SARS viral testing, 797
- Trachoma, 230
- Transbronchial brushing, 586
- Transbronchial lung biopsy, 585
- Transbronchial needle aspiration, 585, 585*f*
- Transesophageal echocardiography (TEE), 340, 903–905, 904*f*
- Transesophageal endoscope, 903, 904*f*
- Transferrin, 546–549, 549*t*  
  normal findings in, 546  
  total iron-binding capacity and, 547
- Transferrin receptor (TfR) assay, 906–907, 907*t*
- Transferrin saturation, 546–549, 549*t*  
  normal findings in, 546  
  total iron-binding capacity and, 547–548
- Transfusion-related acute lung injury (TRALI), 639
- Transfusions, of blood, 153
- Transhepatic cholangiography, percutaneous, 364
- Transthoracic echocardiography (TTE), 339–341
- Transthyretin, 728–729, 729*t*
- Transudates  
  peritoneal, 671  
  in pleural fluid analysis, 864
- Traumatic puncture, 577
- Treponema pallidum*, 859
- Treponema test, 859
- Triamterene  
  urine color and, 932*t*  
  WBC levels and, 976
- Trichomonas vaginalis*, 813*t*, 814  
  in sexually transmitted disease testing, 813–814
- Tricyclics, thromboelastography and, 874–875
- Triglycerides (TGs), 866, 908–909, 908–909*t*
- Triiodothyronine, 910–911, 911*t*
- Triple renal study, 780
- Triple test, 612
- Trisomy 13, 217
- Trisomy 18, 217
- Trisomy 21  
  cell-free maternal DNA testing for, 217  
  maternal screen testing, 612–613  
  maternal serum levels of AFP during, 613
- Troponins, 912–914, 914*t*
- “Trough” anti-Xa tests, 77
- Trough level, drug, 332
- Trypsin, 662–664, 664*t*
- Trypsin-like immunoreactivity, 662
- Trypsinogen, 662
- TSH injection, thyroid-stimulating hormone and, 887
- T-suppressor cells, 219
- Tubal ligation, surgical,  
  hysterosalpingography in, 531
- Tubal obstruction, hysterosalpingography in, 531
- Tuberculin skin test (TST), 915–916
- Tuberculin test, 915–916
- Tuberculosis (TB)  
  blood test, 919  
  culture, 917–918  
  Gold test in, 919–921  
  nucleic acid amplification for, 919–921  
  testing, 919–921, 919*t*  
  tuberculin test for, 915
- Tubular disease, 628
- Tubular (epithelial) casts, 937
- Tularemia, bioterrorism infectious agents testing for, 142*t*, 143
- Tumor analysis  
  breast cancer, 176–178  
  colon cancer, 256–257
- Tumor markers (TMs), 194–196, 195–196*t*  
  in cerebrospinal fluid, 579
- Tumor-associated markers, 194–196, 195–196*t*
- Tumors  
  carcinoid, 807  
  5-HIAA in, 528  
  cardiac, ultrasound for, 339



- Tumors (*Continued*)  
 hormone-producing, 652  
 neuroendocrine, 652  
 rapidly growing, 720–721  
 sputum cytology for, 843  
 stool cancer screening for, 844
- Tumor-specimen analysis, progesterone receptor assay for, 735
- T-wave alternans, microvolt, 344
- Two-dimensional echocardiography, 339
- Type 1 diabetes mellitus (DM), 327
- Type 2 diabetes, insulin assay in, 539
- Type I blocking antibody, 545
- Type I HIT, platelet antibody detection and, 704
- Type I (monoclonal) cryoglobulinemia, 306
- Type I procollagen, amino-terminal propeptide of, 168, 170
- Type II binding antibody, 545
- Type II HIT, platelet antibody detection and, 704
- Tyrosinemia, newborn metabolic screening and, 644
- U**
- U bag, for urine culture and sensitivity, 948
- U wave, 344, 344*f*
- Ulceration, benign peptic, 922
- Ultra-sensitive CRP, 295–296, 296*t*
- Ultrasonography  
 breast, 181–182, 182*t*  
 of scrotum, 799
- Ultrasound  
 abdominal, 1–4, 2*f*  
 of amniotic fluid, 407–408  
 cardiac, 339  
 carotid, 214–215  
 endobronchial, 183–184  
 endourethral urologic, 1  
 fetus and, 685  
 hepatic, 1–2  
 intravascular, 543–544  
 pelvic, 685–687  
 of prostate, 741–742, 741*f*  
 quantitative, in bone mineral density, 157  
 rectal, 741–742, 741*f*  
 scrotal, 799–800, 800*t*  
 tests, 1008  
 thyroid, 891, 891*t*  
 vascular, 960–962
- Ultrasound absorption, in bone mineral density, 157
- Ultrasound mammography, 181–182, 182*t*
- Umbilical artery velocity, 408
- Unconjugated hyperbilirubinemia, 138–139
- Unconjugated (indirect) bilirubin, 137
- Unfractionated standard heparin (UFH), plasma anti-Xa assay in, 77
- Unilateral hearing loss, 357
- Units of measurement, 1015–1016, 1015–1016*t*
- Universal donors, 153
- Universal recipients, 153
- Unsaturated iron-binding capacity (UIBC), 546
- Upper gastrointestinal (UGI) endoscopy, 382–385
- Upper gastrointestinal x-ray study, 922–924
- Upper GI series. *See* Upper gastrointestinal x-ray study
- Upper GI study, air-contrast, 923
- Urea breath test (UBT), 925–926
- Urea crystals, 936
- Urease, 482
- Uremia, 702, 709
- Ureteral catheterization, by cystoscopy, 310*f*
- Ureteral obstruction, 779
- Urethral culture, 814, 816, 817*f*  
 for *Chlamydia*, 231  
 for herpes simplex, 499
- Urethral pressure measurements, 949–952
- Urethral pressure profile (UPP), 949–952
- Uric acid, 927–929, 929*t*
- Uricosuria, 928
- Urinalysis (UA), 930–944, 932*t*, 940–944*t*
- Urinary calcium, 190
- Urinary catheter, indwelling, 948
- Urinary diversion, 948
- Urinary incontinence, pelvic floor sphincter electromyography for, 353
- Urinary stone analysis, 945–946
- Urine  
 amylase in, 56–57  
 appearance of, 931, 937, 940  
 bilirubin in, 935–936, 939  
 bone turnover markers in, 171  
 casts in, 936–937  
 clean-catch collection of, 939, 948  
 collection, in tuberculosis culture, 918  
 color of, 931, 932*t*, 937, 940  
 cortisol in, 290–292, 292*t*  
 crystals in, 936, 939, 943  
 culture, 816  
 delta-aminolevulinic acid, 321, 321*t*  
 11-dTXB2, 710  
 D-xylose in, 983  
 glucose in, 934, 941  
 ketones in, 935, 938–939, 942  
 leucine aminopeptidase in, 561, 561*t*  
 leukocyte esterase in, 935, 938, 942  
 luteinizing hormone in, 593  
 methylmalonic acid in, 969  
 metyrapone test in, 17  
 midstream collection of

## 1062 index

### Urine (*Continued*)

- for urinalysis, 939
  - for urine culture, 948
  - nicotine and metabolites for, 646, 646*t*
  - nitrites in, 935, 942
  - odor of, 931, 937, 940
  - osmolality of, 654, 656–657, 657*t*, 934
  - pepsinogen, 688
  - pH of, 931–933, 937–938, 941
  - platelet function assay and, 709
  - porphyria and, 717
  - potassium, 727, 727*t*
  - prealbumin, 728
  - pregnanediol, 731
  - protein in, 933, 938, 941
  - red blood cells in, 937, 939, 944
  - sodium in, 837–838, 838*t*
  - specific gravity of, 934–935, 938, 942
  - substance abuse testing of, 852
  - suprapubic aspiration of, 948
  - toxicology screening for, 1014
  - uric acid in, 927–929
  - urobilinogen in, 935–936, 939
  - white blood cells in, 937, 939, 944
  - Zika virus in, 985
- Urine culture and sensitivity, 947–948
- Urine drug testing, 852–854, 853*t*
- Urine flow studies, 949–952
- Urine M-spike, monitoring, 134
- Urine osmolar gap, 656
- Urine specimen, for *Chlamydia*, 231
- Urine testing, 392, 1009–1010
- in 17-OCHS, 526
  - for amino acids, 45
  - for HIV, 508
  - for proteins, 134
- Urobilinogen, 935–936, 939
- Urochrome, 931
- Urodynamic studies, 949–952, 949*t*
- Uroflowmetry, 949–952
- Urogenital tract, HSV 2 in, 498
- Urography
- CT, 274
  - intravenous, 764–768
- Urokinase, thrombosis indicators fibrin monomers and, 877
- Urokinase plasminogen activator (uPA), 195–196*t*, 700
- Urologic system, tests for, 997–998
- Urologic tract, ultrasound of, 1
- Uroporphyrinogen-1-synthase, 953, 953*t*
- U.S. Preventive Services Task Force (USPSTF), 668
- mammography screening guidelines, 608–609*t*
- Uterine contractions, 409
- Uterine hyperstimulation, 409
- Uterine problems, hysteroscopy in, 533–534

Uterosalingography, 531–532

Uterotubography, 531–532

## V

- Vagina, ultrasound of, 685–687
- Vaginal secretions, fibronectin and, 412
- Vaginal speculum, 669, 669*f*
- Valproic acid, 333*t*
- Vancomycin, 333*t*
- Vanillylmandelic acid (VMA), 954–957, 957*t*
- Varicella virus testing, 958–959
- Varicella zoster virus (VZV), 958–959
- Varices, barium swallow in, 131
- Vascular ultrasound studies, 960–962
- abnormal findings in, 962
  - interfering factors in, 961
- Vasodilators, 207
- Vasopressin, 71–73, 73*t*
- Vasopressor syncope syndrome, 359
- Velocity, conduction, 355
- Vemurafenib, 257
- Veneral Disease Research Laboratory (VDRL), 859–860
- Venography, of lower extremities, 963–964
- Venous Doppler ultrasound, 960–962
- Venous duplex scan, 960–962
- Venous pressure, central, 202–203*t*
- Ventilation/perfusion scanning, 590–592
- contraindications to, 591
- Ventilation scan, 590–591
- Ventricular natriuretic peptide, 636–637
- Ventricular pressure, left, end-diastolic, 202–203*t*
- Ventriculography, 201–206, 202–203*t*, 204*f*
- Vertebral fracture assessment (VFA), 158
- Vertigo, 357
- Very-low-density lipoproteins (VLDLs), 101–102, 565–568
- apolipoproteins in, 101–102
- Vesicoureteral reflux, 307
- Video colpography, 263–264
- Video colposcopy, 263–264
- Video-assisted thoracotomy (VAT), 870
- Videofluoroscopy swallowing examination, 855–856
- Viral capsid antigen-antibodies (VCAs), 368
- Viral cultures, 965–966, 965*t*
- Viral hepatitis, ALT/AST ratio in, 21
- Viral load
- antiretroviral therapy based on, 504*t*
  - of HBV, 496
  - of HIV, 503–505
  - to predict disease course, 504*t*
- Virtual autopsy, 274
- Virtual bronchoscopy, 281
- Virtual colonoscopy, 273
- Virtual esophagoscopy, 281

Virus testing, 967–968  
 Viscosity, in inflammatory arthritis, 118  
 Visual-evoked responses (VERs), 393–394  
   prolonged latency for, 395  
 Visually read test, 464  
 Vital capacity (VC), 760*f*, 761  
 Vitamin B<sub>6</sub>, dietary deficiency of, in  
   homocysteine, 513  
 Vitamin B<sub>12</sub>, 969–970, 970*t*  
   deficiency of, 426, 969  
   dietary deficiency of, in homocysteine, 513  
 Vitamin D, 694, 971–973, 973*t*  
   calcium and, 189–190  
   deficiency of, 675, 971–972  
   normal findings in, 971  
   procedure and patient care in, 972  
   related physiology in, 971–972  
   test explanation in, 971–972  
   type of test for, 971  
 Vitamin D<sub>2</sub>, 971  
 Vitamin D<sub>3</sub>, 971  
 Vitamin K oxidoreductase (VKOR)  
   complex, 754–755  
 VKOR, prothrombin time, 754–755  
 VKORC1, prothrombin time, 754–755  
 Voiding cystography, 307–308  
 Voiding cystourethrography (VCUG),  
   307–308  
 Volume of isoflow (VisoV), 761–762  
 von Willebrand disease (vWD), 709  
 von Willebrand factor, 248, 249–250*t*  
 V/Q scan, 590–591

## W

Warfarin  
   coumarin derivatives, 754  
   thrombosis indicators fibrin monomers  
     and, 878  
 Warfarin pharmacogenomic test panel,  
   754–755  
 Warfarin sodium, thromboelastography and,  
   874–875  
 Wassermann test, 859  
 Water deprivation test, 71–72  
 Water load test, 72–73  
 Water-soluble contrast material, for  
   myelography, 632  
 Waxy casts, urinary, 936, 943  
 Wedge pressure, pulmonary, 202–203*t*  
 Wegener granulomatosis (WG),  
   antineutrophil cytoplasmic antibody  
     in, 84  
 Western blot assay, for Lyme disease, 595–596  
 Western blot test, 506–509  
 Wetmount preparation, microscopic  
   examination of, of vaginal secretions,  
   814

White blood cells (WBCs)  
   in cerebrospinal fluid, 578  
   lactoferrin and, 555  
   paracentesis and, 671–672  
   in pleural fluid analysis, 865  
   urinary, 937, 939, 944  
 White blood cell (WBC) antigens, 518–519  
 White blood cell (WBC) casts, 937  
 White blood cell (WBC) count, and  
   differential count, 974–978, 974*t*  
   abnormal findings in, 977–978, 977*t*  
   causes of abnormalities in, 977–978*t*  
 White blood cell (WBC) scan, 979–980  
 White blood cell (WBC) stool test, 849  
 White cells, stool for, 849  
 Whole-body thyroid scan, 883  
 Wireless capsule endoscopy, 382  
 Wireless pH probe, 378  
 Wound culture and sensitivity, 981–982

## X

X-ray. *See also* Arteriography (angiography);  
   Computed tomography  
   bone, 172–173, 173*t*  
   in bone densitometry, 159  
   chest, 227–229, 229*t*  
   with contrast dye, 757, 764  
   examinations, 1010  
   KUB (kidney, ureter, bladder), 650  
   lower extremities, 963  
   obstruction series in, 650, 651*t*  
   PET scan and, 719  
   renal pelvis, 764  
   of salivary duct, 818  
   skull, 825  
   small bowel, 830  
   spinal, 839–840  
   supine abdominal study, 650  
   swallowing, 855  
   upper gastrointestinal, 922–924  
 X-ray contrast dye, iodinated, 364  
 Xylose absorption test, 983–984,  
   983–984*t*  
 Xylose tolerance test, 983–984, 983–984*t*

## Y

Y402H, 20  
 Yellow fever, bioterrorism infectious agents  
   testing for, 141, 142*t*  
*Yersinia pestis*, plague in, 142*t*, 143

## Z

Zika virus (ZIKV), 985–986  
 Zinc protoporphyrin (ZPP), 987, 987*t*  
 Zollinger-Ellison (ZE) syndrome, 439  
 Zona pelucida binding tests, 802

## abbreviations for diagnostic and laboratory tests

<b>A</b>	AAT	Alpha <sub>1</sub> -antitrypsin
	ABGs	Arterial blood gases
	ACE	Angiotensin-converting enzyme
	ACT	Activated clotting time
	ACTH	Adrenocorticotrophic hormone
	ADH	Antidiuretic hormone
	AFB	Acid-fast bacilli
	AFP	Alpha-fetoprotein
	A/G ratio	Albumin/globulin ratio
	AIT	Agglutination inhibition test
	ALA	Aminolevulinic acid
	ALP	Alkaline phosphatase
	ALT	Alanine aminotransferase
	AMA	Antimitochondrial antibody
	ANA	Antinuclear antibody
	ANCA	Antineutrophil cytoplasmic antibody
	ANP	Atrial natriuretic peptide
	APCA	Antiparietal cell antibody
	APTT	Activated partial thromboplastin time
	ASMA	Anti-smooth muscle antibody
	ASO	Antistreptolysin O titer
	AST	Aspartate aminotransferase
<b>B</b>	BE	Barium enema
	BMC	Bone mineral content
	BMD	Bone marrow density
	BMP	Basic metabolic panel
	BNP	Brain natriuretic peptide
	BRCA	Breast cancer (gene)
	BSAP	Bone-specific alkaline phosphatase
	BTA	Bladder tumor antigen
	BUN	Blood urea nitrogen
<b>C</b>	C&S	Culture and sensitivity
	CBC	Complete blood count
	CC	Creatinine clearance
	CEA	Carcinoembryonic antigen
	CK	Creatine kinase
	CMP	Comprehensive metabolic panel
	CMV	Cytomegalovirus
	CO	Carbon monoxide
	CO <sub>2</sub>	Carbon dioxide
	COHb	Carboxyhemoglobin test
	CPK, CP	Creatine phosphokinase
	CRP	C-reactive protein
	CSF	Cerebrospinal fluid
	CST	Contraction stress test
	CT	Computed tomography
	cTnI	Cardiac troponin I
	cTnT	Cardiac troponin T
	CVB	Chorionic villus biopsy
	CVS	Chorionic villus sampling
	CXR	Chest x-ray

<b>D</b>	D&C	Dilation and curettage
	DEXA	Dual-energy x-ray absorptiometry
	DSA	Digital subtraction angiography
	DSMA	Disodium monomethane arsonate
	DST	Dexamethasone suppression test
<b>E</b>	ECG, EKG	Electrocardiography
	EEG	Electroencephalogram
	EGD	Esophagogastroduodenoscopy
	EIA	Enzyme immunoassay
	ELISA	Enzyme-linked immunosorbent assay
	EMG	Electromyography
	ENG	Electroneurography
	EP	Evoked potential
	EPCA	Early prostate cancer antigen
	EPO	Erythropoietin
	EPS	Electrophysiologic study
	ER	Estrogen receptor
	ERCP	Endoscopic retrograde cholangiopancreatography
	ESR	Erythrocyte sedimentation rate
	EUG	Excretory urography
<b>F</b>	FBS	Fasting blood sugar
	FDPs	Fibrin degradation products
	%FPSA	Percent free PSA
	FSH	Follicle-stimulating hormone
	FSPs	Fibrin split products
	FT <sub>4</sub>	Thyroxin, free
	FTA-ABS	Fluorescent treponemal antibody absorption test
	FVL	Factor V Leiden
<b>G</b>	GE reflux	Gastroesophageal reflux scan
	GGT	Gamma-glutamyl transferase
	GGTP	Gamma-glutamyl transpeptidase
	GH	Growth hormone
	GHb, GHB	Glycosylated hemoglobin
	GI series	Gastrointestinal series
	GTT	Glucose tolerance test
<b>H</b>	HAA	Hepatitis-associated antigen
	HAI	Hemagglutination inhibition
	Hb, Hgb	Hemoglobin
	HCG	Human chorionic gonadotropin
	HCO <sub>3</sub>	Bicarbonate
	Hct	Hematocrit
	Hcy	Homocysteine
	HDL	High-density lipoprotein
	5-HIAA	Hydroxyindoleacetic acid
	HIDA	Hepatic iminodiacetic acid
	HIV	Human immunodeficiency virus
	HLA-B27	Human lymphocyte antigen B27
	HTLV	Human T-cell lymphotropic virus
<b>I</b>	IAA	Insulin autoantibody
	ICA	Islet cell antibody
	Ig	Immunoglobulin
	INR	International normalized ratio
	IV-GTT	Intravenous glucose tolerance test
	IVP	Intravenous pyelography
	IVU, IUG	Intravenous urography

<b>K</b>	KS	Ketosteroid
	KUB	Kidney, ureter, and bladder
<b>L</b>	LAP	Leucine aminopeptidase
	LATS	Long-acting thyroid stimulator
	LDH	Lactic dehydrogenase
	LDL	Low-density lipoprotein
	LFTs	Liver function tests
	LH	Luteinizing hormone
	LP	Lumbar puncture
	L/S ratio	Lecithin/sphingomyelin ratio
	LS spine	Lumbosacral spine
<b>M</b>	MA	Microalbumin
	MCH	Mean corpuscular hemoglobin
	MCHC	Mean corpuscular hemoglobin concentration
	MCV	Mean corpuscular volume
	MEG	Magnetic encephalography
	M/E ratio	Myeloid/erythroid ratio
	MMA	Methylmalonic acid
	MPG	Mean plasma glucose
	MPV	Mean platelet volume
	MRI	Magnetic resonance imaging
	MUGA	Multigated acquisition cardiac scan
<b>N</b>	NMP 22	Nuclear matrix protein 22
	NST	Nonstress test
	NT <sub>x</sub>	N-telopeptide
<b>O</b>	O&P	Ova and parasites
	OB	Occult blood
	OCT	Oxytocin challenge test
	OGTT	Oral glucose tolerance test
	17-OHCS	17-Hydroxycorticosteroids
<b>P</b>	PAB	Prealbumin
	PAI-1	Plasminogen activator inhibitor-1
	PAP	Prostatic acid phosphatase
	Pco <sub>2</sub>	Partial pressure of carbon dioxide
	PCR	Polymerase chain reaction
	PET	Positron emission tomography
	PFTs	Pulmonary function tests
	pH	Hydrogen ion concentration
	PKU	Phenylketonuria
	PMN	Polymorphonuclear
	PNH	Paroxysmal nocturnal hemoglobinuria
	Po <sub>2</sub>	Partial pressure of oxygen
	PO <sub>4</sub>	Phosphate
	PPBS	Postprandial blood sugar
	PPD	Purified protein derivative
	PPG	Postprandial glucose
	PR	Progesterone receptor
	PRA	Plasma renin assay
	PSA	Prostate-specific antigen
	PT	Prothrombin time
	PTH	Parathormone, parathyroid hormone
	PTHc, PTC	Percutaneous transhepatic cholangiography
	PTT	Partial thromboplastin time
	PYD	Pyridinium crosslink

<b>R</b>	RAST	Radioallergosorbent test
	RBC	Red blood cell
	RDW	Red blood cell distribution width
	RF	Rheumatoid factor
	RIA	Radioimmunoassay
	RPR	Rapid plasma reagin test
<b>S</b>	S&A	Sugar and acetone
	SACE	Serum angiotensin-converting enzyme
	SBF	Small bowel follow-through
	SGOT	Serum glutamic-oxaloacetic transaminase
	SGPT	Serum glutamic-pyruvic transaminase
	SLNB	Sentinel lymph node biopsy
	SPA	Sperm penetration assay
	SPECT	Single-photon emission computed tomography
	STS	Serologic test for syphilis
<b>T</b>	T <sub>3</sub>	Triiodothyronine
	T <sub>4</sub>	Thyroxine
	T&C	Type and crossmatch
	T&S	Type and screen
	TBG	Thyroxine-binding globulin
	TBPA	Thyroxine-binding prealbumin
	TEE	Transesophageal echocardiography
	Tg	Thyroglobulin
	TGs	Triglycerides
	TIBC	Total iron-binding capacity
	TRF	Thyrotropin-releasing factor
	TRH	Thyrotropin-releasing hormone
	TSH	Thyroid-stimulating hormone
	TSI	Thyroid-stimulating immunoglobulins
	TTE	Transthoracic echocardiography
<b>U</b>	UA	Urinalysis
	UGI series	Upper gastrointestinal series
	UPP	Urethral pressure profile
	US	Ultrasound
<b>V</b>	VCUG	Voiding cystourethrography
	VDRL	Venereal Disease Research Laboratory
	VLDL	Very-low-density lipoprotein
	VMA	Vanillylmandelic acid
	VPS	Ventilation/perfusion scanning
<b>W</b>	WBC	White blood cell